

GE Healthcare

# **Optima MR450w**


## **Preinstallation Manual**



OPERATING DOCUMENTATION

5670001  
Revision 2.0

# **WARNING**

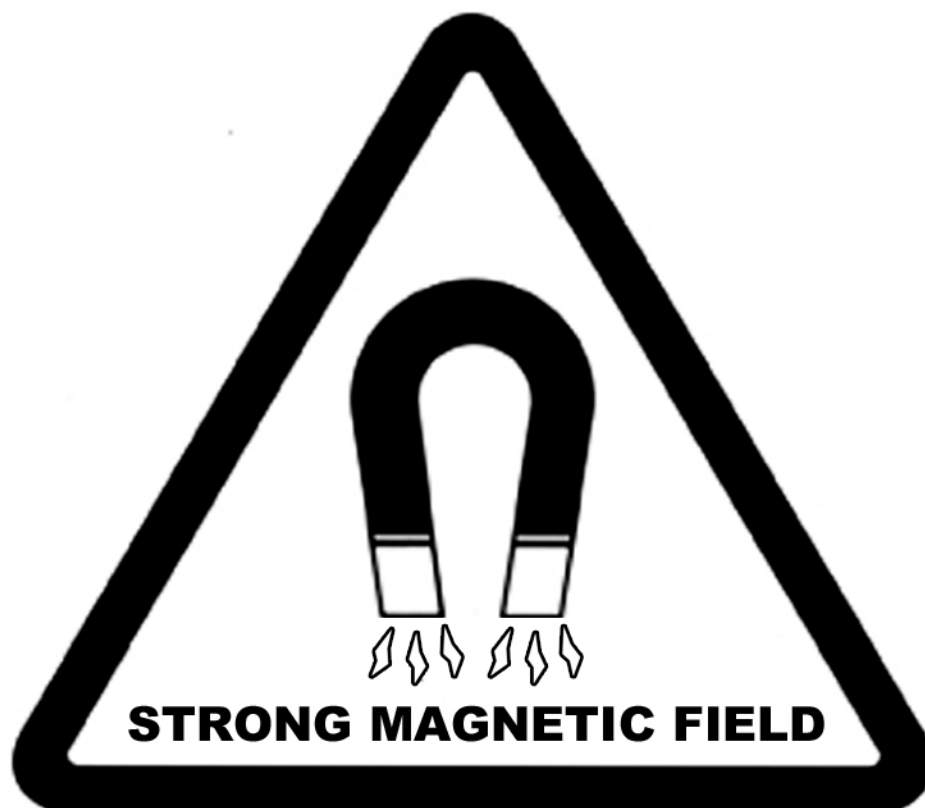

## **NO PACEMAKERS NO METALLIC IMPLANTS**

Persons with pacemakers, neurostimulators or metallic implants must not enter the magnet area. Serious injury may result.



## **NO LOOSE METAL OBJECTS**

Iron and steel materials must not be taken into the magnet area. Serious injury or property damage may result.

# Important Information

## LANGUAGE

### ПРЕДУПРЕЖДЕНИЕ

(BG)

Това упътване за работа е налично само на английски език с изключението на случаите, когато се изисква изрично по друг начин от местното законодателство, или от споразумение на местно ниво.

- Ако доставчикът на услугата на клиента изиска друг език, задължение на клиента е да осигури превод.
- Не използвайте оборудването, преди да сте се консултирали и разбрали упътването за работа.
- Неспазването на това предупреждение може да доведе до нараняване на доставчика на услугата, оператора или пациента в резултат на токов удар, механична или друга опасност.

### 警告

(ZH-CN)

本维修手册仅提供英文版本，除非当地法律或本地供应协议另有明确要求。

- 如果客户的维修服务人员需要非英文版本，则客户需自行提供翻译服务。
- 未详细阅读和完全理解本维修手册之前，不得进行维修。
- 忽略本警告可能对维修服务人员、操作人员或患者造成触电、机械伤害或其他形式的伤害。

### VÝSTRAHA

(CS)

Tento provozní návod existuje pouze v anglickém jazyce, pokud není jinak výslovně vyžadováno místními zákony nebo odsouhlaseno na místní úrovni.

- V případě, že externí služba zákazníkům potřebuje návod v jiném jazyce, je zajištění překladu do odpovídajícího jazyka úkolem zákazníka.
- Nesnažte se o údržbu tohoto zařízení, aniž byste si přečetli tento provozní návod a pochopili jeho obsah.
- V případě nedodržování této výstrahy může dojít k poranění pracovníka prodejního servisu, obslužného personálu nebo pacientů vlivem elektrického proudu, respektive vlivem mechanických či jiných rizik.

**ADVARSEL**

(DA)

Denne servicemanual findes kun på engelsk, medmindre andet kræves i henhold til lokal lovgivning eller lokal aftale.

- Hvis en kundes tekniker har brug for et andet sprog end engelsk, er det kundens ansvar at sørge for oversættelse.
- Forsøg ikke at servicere udstyret uden at læse og forstå denne servicemanual.
- Manglende overholdelse af denne advarsel kan medføre skade på grund af elektrisk stød, mekanisk eller anden fare for teknikeren, operatøren eller patienten.

**WAARSCHUWING**

(NL)

Deze onderhoudshandleiding is enkel in het Engels verkrijgbaar, tenzij expliciet vereist door plaatselijke regelgeving of overeengekomen op lokaal niveau.

- Als het onderhoudspersoneel een andere taal vereist, dan is de klant verantwoordelijk voor de vertaling ervan.
- Probeer de apparatuur niet te onderhouden alvorens deze onderhoudshandleiding werd geraadpleegd en begrepen is.
- Indien deze waarschuwing niet wordt opgevolgd, zou het onderhoudspersoneel, de operator of een patiënt gewond kunnen raken als gevolg van een elektrische schok, mechanische of andere gevaren.

**WARNING**

(EN)

This service manual is available in english only except as otherwise expressly required by local law or agreed to at a local level.

- If a customer's service provider requires a language other than english, it is the customer's responsibility to provide translation services.
- Do not attempt to service the equipment unless this service manual has been consulted and is understood.
- Failure to heed this warning may result in injury to the service provider, operator or patient from electric shock, mechanical or other hazards.



**HOIATUS**

(ET)

See teenindusjuhend on saadaval ainult inglise keeles, kui kohalikud seadused ei ütle teisiti või kui kohalikes õigusaktides ei ole otseselt teisiti ette nähtud.

- Kui klienditeeninduse osutaja nõuab juhendit inglise keelest erinevas keeles, vastutab klient tõlketeenuse osutamise eest.
- Ärge üritage seadmeid teenindada enne eelnevalt käesoleva teenindusjuhendiga tutvumist ja sellest aru saamist.
- Käesoleva hoiatuse eiramine võib põhjustada teenuseosutaja, operaatori või patsiendi vigastamist elektrilöögi, mehaanilise või muu ohu tagajärjel.

**VAROITUS**

(FI)

Tämä huolto-ohje on saatavilla vain englanniksi, ellei paikallinen laki nimenomaan toisin vaadi tai jos toisin on sovittu paikallisella tasolla.

- Jos asiakkaan huoltohenkilöstö vaatii muuta kuin englanninkielistä materiaalia, tarvittavan käännöksen hankkiminen on asiakkaan vastuulla.
- Älä yritä korjata laitteistoa ennen kuin olet varmasti lukenut ja ymmärtänyt tämän huolto-ohjeen.
- Mikäli tätä varoitusta ei noudateta, seurauksena voi olla huoltohenkilöstön, laitteiston käyttäjän tai potilaan vahingoittuminen sähköiskun, mekaanisen vian tai muun vaaratilanteen vuoksi.

**ATTENTION**

(FR)

Sauf exigence contraire des lois locales ou accord contraire au niveau local, ce manuel d'installation et de maintenance n'est disponible qu'en anglais.

- Si le technicien d'un client a besoin de ce manuel dans une langue autre que l'anglais, il incombe au client de le faire traduire.
- Ne pas tenter d'intervenir sur les équipements tant que ce manuel d'installation et de maintenance n'a pas été consulté et compris.
- Le non-respect de cet avertissement peut entraîner chez le technicien, l'opérateur ou le patient des blessures dues à des dangers électriques, mécaniques ou autres.

**WARNUNG**

(DE)

Diese Serviceanleitung existiert nur in englischer Sprache, sofern nichts anderes gesetzlich vorgeschrieben oder auf lokaler Ebene vereinbart wurde.

- Falls ein fremder Kundendienst eine andere Sprache benötigt, ist es Aufgabe des Kunden für eine Entsprechende Übersetzung zu sorgen.
- Versuchen Sie nicht diese Anlage zu warten, ohne diese Serviceanleitung gelesen und verstanden zu haben.
- Wird diese Warnung nicht beachtet, so kann es zu Verletzungen des Kundendiensttechnikers, des Bedieners oder des Patienten durch Stromschläge, mechanische oder sonstige Gefahren kommen.

**ΠΡΟΕΙΔΟΠΟΙΗΣΗ**

(EL)

Το παρόν εγχειρίδιο σέρβις διατίθεται μόνο στα αγγλικά, εκτός αν η τοπική νομοθεσία απαιτεί κάτι άλλο ή υπάρχει διαφορετική συμφωνία σε τοπικό επίπεδο.

- Εάν ο τεχνικός σέρβις ενός πελάτη απαιτεί το παρόν εγχειρίδιο σε γλώσσα εκτός των αγγλικών, αποτελεί ευθύνη του πελάτη να παρέχει τις υπηρεσίες μετάφρασης.
- Μην επιχειρήσετε την εκτέλεση εργασιών σέρβις στον εξοπλισμό αν δεν έχετε συμβουλευτεί και κατανοήσει το παρόν εγχειρίδιο σέρβις.
- Αν δεν προσέξετε την προειδοποίηση αυτή, ενδέχεται να προκληθεί τραυματισμός στον τεχνικό σέρβις, στο χειριστή ή στον ασθενή από ηλεκτροπληξία, μηχανικούς ή άλλους κινδύνους.

**FIGYELMEZTETÉS**

(HU)

Ezen karbantartási kézikönyv kizárólag angol nyelven érhető el, kivéve ha a helyi rendelkezések ezt kifejezetten elő nem írják, illetve ha helyi szinten erről külön megállapodás nem születik.

- Ha a vevő szolgáltatója angoltól eltérő nyelvre tart igényt, akkor a vevő felelőssége a fordítás elkészítése.
- Ne próbálja elkezdni használni a berendezést, amíg a karbantartási kézikönyvben leírtakat nem értelmezték.
- Ezen figyelmeztetés figyelmen kívül hagyása a szolgáltató, működtető vagy a beteg áramütés, mechanikai vagy egyéb veszélyhelyzet miatti sérülését eredményezheti.

**AÐVÖRUN**

(IS)

Þessi þjónustuhandbók er eingöngu fánleg á ensku, nema annað sé sérstaklega krafist, löglega eða samþykkt á landsgrundvelli.

- Ef að þjónustuveitandi viðskiptamanns þarfnast annas tungumáls en ensku, er það skylda viðskiptamanns að skaffa tungumálaþjónustu.
- Reynið ekki að afgreiða tækið nema að þessi þjónustuhandbók hefur verið skoðuð og skilin.
- Brot á sinna þessari aðvörun getur leitt til meiðsla á þjónustuveitanda, stjórnanda eða sjúklings frá raflosti, vélrænu eða öðrum áhættum.

**AVVERTENZA**

(IT)

Il presente manuale di manutenzione è disponibile soltanto in inglese, eccetto quando espressamente richiesto dalle normative locali o convenuto a livello locale.

- Se un addetto alla manutenzione richiede il manuale in una lingua diversa, il cliente è tenuto a provvedere direttamente alla traduzione.
- Procedere alla manutenzione dell'apparecchiatura solo dopo aver consultato il presente manuale ed averne compreso il contenuto.
- Il mancato rispetto della presente avvertenza potrebbe causare lesioni all'addetto alla manutenzione, all'operatore o ai pazienti provocate da scosse elettriche, urti meccanici o altri rischi.

**警告**

(JA)

このサービスマニュアルには英語版しかありません。ただし使用国の法令に別異の定めがある、あるいは現地で別段の合意がある場合を除きます。

- サービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないでください。
- この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

**경고**

(KO)

현지 법률에 따라 명시적으로 요구하거나 현지 수준에서 합의한 경우를 제외하고 본 서비스 매뉴얼은 영어로만 이용하실 수 있습니다.

- 고객의 서비스 제공자가 영어 이외의 언어를 요구할 경우, 번역 서비스를 제공하는 것은 고객의 책임입니다.
- 본 서비스 매뉴얼을 참조하여 숙지하지 않은 이상 해당 장비를 수리하려고 시도하지 마십시오.
- 본 경고 사항에 유의하지 않으면 전기 쇼크, 기계적 위험, 또는 기타 위험으로 인해 서비스 제공자, 사용자 또는 환자에게 부상을 입힐 수 있습니다.

**BRĪDINĀJUMS**

(LV)

Šī apkalpes rokasgrāmata ir pieejama tikai angļu valodā, izņemot gadījumus, kad vietējie likumi nepārprotami nosaka citādi vai panākta vienošanās vietējā līmenī.

- Ja klienta apkalpes sniedzējam nepieciešama informācija citā valodā, nevis angļu, klienta pienākums ir nodrošināt tulkošanu.
- Neveiciet aprīkojuma apkalpi bez apkalpes rokasgrāmatas izlasīšanas un saprašanas.
- Šī brīdinājuma neievērošana var radīt elektriskās strāvas trieciena, mehānisku vai citu risku izraisītu traumu apkalpes sniedzējam, operatoram vai pacientam.

**ĮSPĖJIMAS**

(LT)

Šis eksploatavimo vadovas yra tik anglų kalba, išskyrus tuos atvejus, kai vietiniai įstatymai tiesiogiai numato kitokius reikalavimus arba vietinių lygiu sutarta kitaip.

- Jei kliento paslaugų tiekėjas reikalauja vadovo kita kalba – ne anglų, suteikti vertimo paslaugas privalo klientas.
- Nemėginkite atlikti įrangos techninės priežiūros, jei neperskaitėte ar nesupratote šio eksploatavimo vadovo.
- Jei nepaisysite šio įspėjimo, galimi paslaugų tiekėjo, operatoriaus ar paciento sužalojimai dėl elektros šoko, mechaninių ar kitų pavojų.

**ADVARSEL**

(NO)

Denne servicehåndboken finnes bare på engelsk, bortsett fra dersom det motsatte uttrykkelig er fastsatt av lokal lovgivning eller det er inngått annen avtale lokalt.

- Hvis kundens serviceleverandør trenger et annet språk, er det kundens ansvar å sørge for oversettelse.
- Ikke forsøk å reparere utstyret uten at denne servicehåndboken er lest og forstått.
- Manglende hensyn til denne advarselen kan føre til at serviceleverandøren, operatøren eller pasienten skades på grunn av elektrisk støt, mekaniske eller andre farer.

**OSTRZEŻENIE**

(PL)

Niniejszy podręcznik serwisowy dostępny jest jedynie w języku angielskim, chyba że lokalne przepisy lub umowy wyraźnie stanowią inaczej.

- Jeśli dostawca usług klienta wymaga języka innego niż angielski, zapewnienie usługi tłumaczenia jest obowiązkiem klienta.
- Nie próbować serwisować wyposażenia bez zapoznania się z niniejszym podręcznikiem serwisowym i zrozumienia go.
- Niezastosowanie się do tego ostrzeżenia może spowodować urazy dostawcy usług, operatora lub pacjenta w wyniku porażenia prądem elektrycznym, zagrożenia mechanicznego bądź innego.

**ATENȚIE**

(RO)

Acest manual de service este disponibil numai în limba engleză, cu excepția cazului în care este o cerință obligatorie stipulată de legislația națională sau convenită la nivel local.

- Dacă un furnizor de servicii pentru clienți necesită o altă limbă decât cea engleză, este de datoria clientului să furnizeze o traducere.
- Nu încercați să reparați echipamentul decât ulterior consultării și înțelegerii acestui manual de service.
- Ignorarea acestui avertisment ar putea duce la rănirea depanatorului, operatorului sau pacientului în urma pericolelor de electrocutare, mecanice sau de altă natură.

**ОСТОРОЖНО!**

(RU)

Данное руководство по техническому обслуживанию предлагается только на английском языке, за исключением тех случаев, когда наличие руководства на национальном языке является требованием местного законодательства или когда выпуск такого руководства согласован с местным представительством.

- Если сервисному персоналу клиента необходимо руководство не на английском, а на каком-то другом языке, клиенту следует самостоятельно обеспечить перевод.
- Перед техническим обслуживанием оборудования обязательно обратитесь к данному руководству и поймите изложенные в нем сведения.
- Несоблюдение требований данного предупреждения может привести к тому, что специалист по техобслуживанию, оператор или пациент получит удар электрическим током, механическую травму или другое повреждение.

**UPOZORNENIE**

(SK)

Tento návod na obsluhu je k dispozícii len v angličtine, okrem prípadov, kedy tak výslovne vyžadujú miestne zákony alebo je dohodnuté na miestnej úrovni.

- Ak zákaznikov poskytovateľ služieb vyžaduje iný jazyk ako angličtinu, poskytnutie prekladateľských služieb je zodpovednosťou zákazníka.
- Nepokúšajte sa o obsluhu zariadenia, kým si neprečítate návod na obsluhu a neporozumiete mu.
- Zanedbanie tohto upozornenia môže spôsobiť zranenie poskytovateľa služieb, obsluhujúcej osoby alebo pacienta elektrickým prúdom, mechanické alebo iné ohrozenie.

**ATENCION**

(ES)

Este manual de servicio sólo existe en inglés, salvo que la legislación local exija de forma expresa lo contrario, o así se haya acordado a nivel local.

- Si el encargado de mantenimiento de un cliente necesita un idioma que no sea el inglés, el cliente deberá encargarse de la traducción del manual.
- No se deberá dar servicio técnico al equipo, sin haber consultado y comprendido este manual de servicio.
- La no observancia del presente aviso puede dar lugar a que el proveedor de servicios, el operador o el paciente sufran lesiones provocadas por causas eléctricas, mecánicas o de otra naturaleza.

**VARNING**

(SV)

Den här servicehandboken finns bara tillgänglig på engelska om inte annat uttryckligen krävs av lokal lag eller har överenskommit på lokal nivå.

- Om en kunds servicetekniker har behov av ett annat språk än engelska, ansvarar kunden för att tillhandahålla översättningstjänster.
- Försök inte utföra service på utrustningen om du inte har läst och förstår den här servicehandboken.
- Om du inte tar hänsyn till den här varningen kan det resultera i skador på serviceteknikern, operatören eller patienten till följd av elektriska stötar, mekaniska faror eller andra faror.

**DİKKAT**

(TR)

Aksi, yerel bir yasa tarafından açıkça gerekli görülmediği veya yerel bir seviyede kabul edilmediği takdirde, bu servis kılavuzunun sadece İngilizcesi mevcuttur.

- Eğer müşteri teknisyeni bu kılavuzu İngilizce dışında bir başka lisandan talep ederse, bunu tercüme ettirmek müşteriye düşer.
- Servis kılavuzunu okuyup anlamadan ekipmanlara müdahale etmeyiniz.
- Bu uyarıya uyulmaması, elektrik, mekanik veya diğer tehlikelerden dolayı teknisyen, operatör veya hastanın yaralanmasına yol açabilir.

**AVISO**

(PT-BR)

Este manual de assistência técnica encontra-se disponível unicamente em inglês, salvo disposições em contrário previstas pela legislação local ou acordadas no âmbito local.

- Se outro serviço de assistência técnica solicitar a tradução deste manual, caberá ao cliente fornecer os serviços de tradução.
- Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica.
- A não observância deste aviso pode ocasionar ferimentos no técnico, operador ou paciente decorrentes de choques elétricos, mecânicos ou outros.

**ATENÇÃO**  
(PT-PT)

Este manual de assistência técnica só se encontra disponível em inglês, salvo requisição expressa pela legislação local ou acordo efectuado a nível local.

- Se qualquer outro serviço de assistência técnica solicitar este manual noutro idioma, é da responsabilidade do cliente fornecer os serviços de tradução.
- Não tente reparar o equipamento sem ter consultado e compreendido este manual de assistência técnica.
- O não cumprimento deste aviso pode colocar em perigo a segurança do técnico, do operador ou do paciente devido a choques eléctricos, mecânicos ou outros.

**UPOZORENJE**  
(SR)

Ovo servisno uputstvo je dostupno samo na engleskom jeziku, sem ako lokalni zakon to izričito zahteva ili je dogovoreno na lokalnom nivou.

- Ako klijentov serviser zahteva neki drugi jezik, klijent je dužan da obezbedi prevodilačke usluge.
- Ne pokušavajte da opravite uređaj ako niste pročitali i razumeli ovo servisno uputstvo.
- Zanemarivanje ovog upozorenja može dovesti do povređivanja serviser, rukovaoca ili pacijenta usled strujnog udara ili mehaničkih i drugih opasnosti.

**UPOZORENJE**  
(HR)

Ovaj servisni priručnik dostupan je na engleskom jeziku.

- Ako davatelj usluge klijenta treba neki drugi jezik, klijent je dužan osigurati prijevod.
- Ne pokušavajte servisirati opremu ako niste u potpunosti pročitali i razumjeli ovaj servisni priručnik.
- Zanemarite li ovo upozorenje, može doći do ozljede davatelja usluge, operatera ili pacijenta uslijed strujnog udara, mehaničkih ili drugih rizika.

**警告**  
(ZH-TW)

本維修手冊僅有英文版。

- 若客戶的維修廠商需要英文版以外的語言，應由客戶自行提供翻譯服務。
- 請勿試圖維修本設備，除非 您已查閱並瞭解本維修手冊。
- 若未留意本警告，可能導致維修廠商、操作員或病患因觸電、機械或其他危險而受傷。

**警告**  
(ZH-HK)

本服務手冊僅提供英文版本。

- 倘若客戶的服務供應商需要英文以外之服務手冊，客戶有責任提供翻譯服務。
- 除非已參閱本服務手冊及明白其內容，否則切勿嘗試維修設備。
- 不遵從本警告或會令服務供應商、網絡供應商或病人受到觸電、機械性或其他危險。



## Revision History

| Revision | Date       | Description                                       |                                                               |                                                                                                                                              |
|----------|------------|---------------------------------------------------|---------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------------|
| 1.0      | 01July2009 | Initial Release                                   |                                                               |                                                                                                                                              |
| 2.0      | 01Dec2009  | Document review and approval per DOC0675763 R2 V3 |                                                               |                                                                                                                                              |
|          |            | <b>iTrak No.</b>                                  | <b>Ch, Sec No.</b>                                            | <b>Change</b>                                                                                                                                |
|          |            | 13252239                                          | Ch 6, Sec 1.3                                                 | Added Brainwave cable information                                                                                                            |
|          |            | 13230415                                          | Ch 3, Sec 8.2                                                 | Clarified Incoming air must contain 5% outside air (from outside the room)                                                                   |
|          |            | 13254391                                          | Ch 4, Sec 5                                                   | Updated Pen cabinet installed height and width                                                                                               |
|          |            | <b>SPR No.</b>                                    |                                                               |                                                                                                                                              |
|          |            | MRlhc45228;<br>MRlhc45709                         | Ch 6, Sec 3                                                   | Updated system cable list                                                                                                                    |
|          |            | MRlhc45702                                        | Ch 2, Sec 11                                                  | Standby (no scan) missing, should be 17 kVA                                                                                                  |
|          |            |                                                   | Ch 2, Sec 3                                                   | Missing dimension on illustration 2-2 (Front of Magnet to front of service area = 103.5 in. )                                                |
|          |            |                                                   | Ch 2, Sec 8.2                                                 | 79.2 in. is shown as is 3000 mm. Should be 79 (2000)                                                                                         |
|          |            |                                                   | Ch 3, Sec 5.4                                                 | Missing RF survey specifications (same as 450). Add back in.                                                                                 |
|          |            |                                                   | Ch 3, Sec 5.6                                                 | Magnet anchor drawing titled "Magnet and Dock Anchor Mounting Details"                                                                       |
|          |            | MRlhc45703                                        | Ch 3, Sec 2.2                                                 | Updated steel limit table (distance from isocenter values)                                                                                   |
|          |            | MRlhc45705                                        | Ch 3, Sec 5.2.5                                               | Updated RF shielding test requirement from 150 MHz to 100 MHz (with recommendation for 150 MHz for new construction to accommodate upgrades) |
|          |            | MRlhc45706                                        | Ch 4, Sec 4; Ch 4, Sec 5                                      | Updated gauss limits for PEN and HEC cabinets to 50 Gauss                                                                                    |
|          |            | MRlhc45707                                        | Ch 2, Sec 4                                                   | Updated EMC declaration per 60601-1-2 Edition 2.1                                                                                            |
|          |            | MRlhc45708                                        | Ch 3, Sec 7.2                                                 | Updated Table Illustration to show DV table (minor dimension changes)                                                                        |
|          |            | MRlhc45709                                        | Ch 2, Sec 9; Ch 6, Sec 4                                      | Added/updated Facility cooling requirements and hose/hose insulation specifications                                                          |
|          |            | MRlhc45714                                        | Ch 3, Sec 9.3                                                 | Updated Common Ground Stud requirement to increase allowable filter and pipe distances                                                       |
|          |            | MRlhc45715                                        | Ch 2, Sec 11                                                  | Changed MDP Regulation specification to 2%                                                                                                   |
|          |            | MRlhc45718                                        | Ch 3, Sec 6.4                                                 | Added Vibromat dimensions to allow floor loading calculations                                                                                |
|          |            | MRlhc46104                                        | Ch 2, Sec 12.3                                                | PEN, HEC, PGR, SPW cabinet shipping dimensions updated                                                                                       |
|          |            | MRlhc46124                                        | Ch 3, Sec 6.3.1                                               | Added ceiling interface drawing and requirements                                                                                             |
|          |            | <b>Other</b>                                      |                                                               |                                                                                                                                              |
|          |            |                                                   | Ch 2, Sec 10                                                  | Illustration 2-15, changed "Top of HEC" to "Terminal Strip" and remove E3027                                                                 |
|          |            | ECO 2076814                                       | Ch 3, Sec 7.3                                                 | New MRU specifications and requirements added                                                                                                |
|          |            |                                                   | Ch 4, Sec 2                                                   | New MDP Updates                                                                                                                              |
|          |            |                                                   | Ch 2, Sec 2.2; Ch 3, Sec 5.2.6; Ch 4; Ch 6 Sec 2; Ch 6, Sec 4 | Added MRE requirements and specifications                                                                                                    |
|          |            |                                                   | Ch 4, Sec 7                                                   | Updated PEN Panel mounting drawing                                                                                                           |

This page left intentionally blank.

# Table of Contents

|                                                                 |           |
|-----------------------------------------------------------------|-----------|
| <b>Chapter 1 INTRODUCTION.....</b>                              | <b>21</b> |
| 1 Preinstall Manual Introduction.....                           | 21        |
| 1.1 Document Purpose.....                                       | 21        |
| 1.2 Intended User.....                                          | 21        |
| 1.3 Document Overview.....                                      | 21        |
| <b>Chapter 2 GENERAL SYSTEM LEVEL.....</b>                      | <b>23</b> |
| 1 Upgrade Requirements.....                                     | 23        |
| 2 System Components.....                                        | 24        |
| 2.1 Magnet Room.....                                            | 24        |
| 2.2 Equipment Room.....                                         | 24        |
| 2.3 Control Room.....                                           | 24        |
| 2.4 Accessories.....                                            | 24        |
| 3 MR Suite Minimum Room Size Requirements.....                  | 25        |
| 4 IEC EMC Compliance.....                                       | 28        |
| 5 MR System Seismic Requirements.....                           | 29        |
| 6 MR Suite Acoustic Specifications.....                         | 30        |
| 6.1 Acoustic Specifications.....                                | 30        |
| 6.2 Structureborne Vibration Control Specifications.....        | 30        |
| 7 MR Suite Magnetic Field Specifications.....                   | 31        |
| 7.1 Magnetic Fringe Field.....                                  | 31        |
| 7.2 Interference from Changing Magnetic Fields.....             | 33        |
| 7.3 Electrical Current.....                                     | 36        |
| 7.4 Non-MR System Equipment Sensitivity to Magnetic Fields..... | 37        |
| 8 Multiple MR System Requirements.....                          | 39        |
| 8.1 Multiple Magnets.....                                       | 39        |
| 8.2 Shared Equipment Rooms.....                                 | 39        |
| 9 MR Suite Temperature and Humidity.....                        | 41        |
| 9.1 Temperature and Humidity Requirements.....                  | 41        |
| 9.2 Equipment Heat Output Specifications.....                   | 41        |

|                                                                               |           |
|-------------------------------------------------------------------------------|-----------|
| 10 Facility Coolant Requirements.....                                         | 43        |
| 10.1 Heat Exchanger Cabinet (HEC) Coolant Requirements.....                   | 43        |
| 10.2 Cryocooler Compressor (CRY) Emergency Facility Coolant Requirements..... | 44        |
| 11 MR Suite Electrical Requirements.....                                      | 47        |
| 12 MR System Shipping and Receiving.....                                      | 50        |
| 12.1 MR System Shipping Requirements.....                                     | 50        |
| 12.2 Facility Delivery Route Requirements.....                                | 50        |
| 12.3 MR System Component Shipping Specifications.....                         | 50        |
| <b>Chapter 3 MAGNET ROOM.....</b>                                             | <b>53</b> |
| 1 Magnet Room Introduction.....                                               | 53        |
| 2 Magnet Room Structural Requirements.....                                    | 54        |
| 2.1 Overview.....                                                             | 54        |
| 2.2 Environmental Steel Limits.....                                           | 54        |
| 2.3 Vibration Requirements.....                                               | 55        |
| 3 Magnetic Shielded Room Requirements.....                                    | 57        |
| 4 Acoustic Room Specifications.....                                           | 58        |
| 5 RF Shielded Room Requirements.....                                          | 59        |
| 5.1 RF Noise Definitions.....                                                 | 59        |
| 5.2 RF Shield Requirements.....                                               | 59        |
| 5.2.1 Customer Responsibilities.....                                          | 59        |
| 5.2.2 Requirements.....                                                       | 59        |
| 5.3 Broadband Interference Avoidance Requirements.....                        | 60        |
| 5.4 Ambient Radio Frequency Interference (RFI) Requirements.....              | 60        |
| 5.5 Magnet Mounting Requirements.....                                         | 61        |
| 5.6 Dock Anchor Mounting Requirements.....                                    | 62        |
| 6 Finished Room Requirements.....                                             | 64        |
| 6.1 Walls.....                                                                | 64        |
| 6.2 Doors and Magnet Access Openings.....                                     | 64        |
| 6.3 Finished Ceiling.....                                                     | 65        |
| 6.4 Magnet Room Floors.....                                                   | 65        |

|                  |                                                                |           |
|------------------|----------------------------------------------------------------|-----------|
| 7                | Magnet Room Equipment Specifications.....                      | 67        |
| 7.1              | Magnet (MAG) Specifications.....                               | 67        |
| 7.2              | Patient Table (PT) Specifications.....                         | 67        |
| 7.3              | Magnet Rundown Unit (MRU) Specifications and Requirements..... | 68        |
| 7.4              | Oxygen Monitor Sensor Specifications.....                      | 70        |
| 8                | Magnet Room Venting Requirements.....                          | 71        |
| 8.1              | Overall Requirements.....                                      | 71        |
| 8.2              | HVAC Vent Requirements.....                                    | 71        |
| 8.3              | Emergency Exhaust Vent Requirements.....                       | 71        |
| 8.4              | Pressure Equalization Vent Requirement.....                    | 73        |
| 8.5              | Cryogenic Venting.....                                         | 73        |
| 8.5.1            | Interior Vent Requirements.....                                | 73        |
| 8.5.1.1          | General.....                                                   | 73        |
| 8.5.1.2          | Vent Size.....                                                 | 73        |
| 8.5.1.3          | Vent Materials.....                                            | 74        |
| 8.5.1.4          | Support.....                                                   | 74        |
| 8.5.1.5          | Construction.....                                              | 74        |
| 8.5.2            | Exterior Vent Requirements.....                                | 75        |
| 8.5.2.1          | Vent Support.....                                              | 75        |
| 8.5.2.2          | Vent Construction.....                                         | 75        |
| 8.5.2.3          | Vent Exit.....                                                 | 75        |
| 9                | Magnet Room Electrical and Grounding Requirements.....         | 79        |
| 9.1              | Electrical Line and Filter Requirements.....                   | 79        |
| 9.2              | Lighting Requirements.....                                     | 79        |
| 9.3              | Grounding Requirements.....                                    | 79        |
| <b>Chapter 4</b> | <b>EQUIPMENT ROOM.....</b>                                     | <b>81</b> |
| 1                | Equipment Room Overview.....                                   | 81        |
| 2                | Main Disconnect Panel (MDP) Specifications.....                | 82        |
| 3                | Power, Gradient, RF Cabinet (PGR) Specifications.....          | 83        |
| 4                | Heat Exchanger Cabinet (HEC) Specifications.....               | 84        |

|                                                              |            |
|--------------------------------------------------------------|------------|
| 5 Penetration Panel Cabinet (PEN) Specifications.....        | 86         |
| 6 Secondary Penetration Wall (SPW) Specifications.....       | 88         |
| 7 PEN and SPW Wall Opening Requirements.....                 | 89         |
| 8 Cryocooler Compressor (CRY) Specifications.....            | 90         |
| 9 Magnet Monitor (MON) Requirements and Specifications.....  | 91         |
| 9.1 Requirements.....                                        | 91         |
| 9.2 Specifications.....                                      | 91         |
| 10 Brainwave Lite (BW) Specifications.....                   | 93         |
| 11 Magnetic Resonance Elastography (MRE) Specifications..... | 94         |
| 11.1 Requirements.....                                       | 94         |
| 11.2 Specifications.....                                     | 94         |
| <b>Chapter 5 CONTROL ROOM.....</b>                           | <b>95</b>  |
| 1 Operator Workspace Equipment Specifications.....           | 95         |
| 1.1 Operator Workspace (OW).....                             | 95         |
| 1.2 Global Operator Cabinet (GOC).....                       | 95         |
| 1.3 Host Display.....                                        | 96         |
| 1.4 Host Keyboard.....                                       | 97         |
| 1.5 Pneumatic Patient Alert.....                             | 97         |
| 2 Oxygen Monitor (OXY) Specifications.....                   | 99         |
| 2.1 Oxygen Monitor (OXY) Requirements.....                   | 99         |
| 2.2 Oxygen Monitor Sensor.....                               | 101        |
| <b>Chapter 6 MR SYSTEM INTERCONNECTS.....</b>                | <b>103</b> |
| 1 MR System Interconnects Specifications.....                | 103        |
| 1.1 Component Designator Definitions.....                    | 103        |
| 1.2 Usable Cable Lengths.....                                | 103        |
| 1.3 Brainwave Option.....                                    | 104        |
| 1.4 Magnetic Resonance Elastography (MRE) Option.....        | 104        |
| 2 MR System Interconnects Routing Requirements.....          | 106        |
| 2.1 General Requirements.....                                | 106        |
| 2.2 Magnet Room Requirements.....                            | 106        |

|                                                              |            |
|--------------------------------------------------------------|------------|
| 2.3 Equipment Room Requirements.....                         | 109        |
| 3 MR System Cable Specifications.....                        | 110        |
| 4 Facility Supplied System Interconnects Specifications..... | 119        |
| <b>Chapter 7 APPENDICIES.....</b>                            | <b>121</b> |
| 1 Glossary.....                                              | 121        |
| 2 MR Site Vibration Test Guidelines.....                     | 124        |
| 2.1 Test Measurements.....                                   | 124        |
| 2.2 Equipment (Spectral Analyzer) Set-Up.....                | 124        |
| 2.3 Data Collection.....                                     | 124        |
| 2.3.1 Ambient Baseline Condition.....                        | 124        |
| 2.3.2 Normal Condition.....                                  | 125        |
| 2.4 Presentation/Interpretation of Results.....              | 125        |
| 3 RF Shielding Background and Guidelines.....                | 131        |
| 3.1 RF Shielding Background.....                             | 131        |
| 3.1.1 Discrete Frequency.....                                | 131        |
| 3.1.2 Broadband RF Noise.....                                | 131        |
| 3.2 RF Shielded Room Grounding.....                          | 131        |
| 4 RF Shielded Enclosure Test Guidelines.....                 | 132        |
| 4.1 Introduction.....                                        | 132        |
| 4.1.1 Purpose of Test Plan.....                              | 132        |
| 4.1.2 Not used.....                                          | 132        |
| 4.2 Applicable Documents.....                                | 132        |
| 4.3 Test Sample Set-Up.....                                  | 132        |
| 4.4 Shielding Effectiveness.....                             | 132        |
| 4.5 Measurement Procedure.....                               | 133        |
| 4.5.1 Test Position.....                                     | 133        |
| 4.5.2 Frequency Range.....                                   | 133        |
| 4.5.3 Free Field Calibration.....                            | 133        |
| 4.6 Enclosure Power Reference Isolation.....                 | 133        |
| 4.7 Test Equipment.....                                      | 134        |

|                                                                |     |
|----------------------------------------------------------------|-----|
| 4.8 Data Recording and Verification.....                       | 134 |
| 4.9 Test Report.....                                           | 134 |
| 5 Acoustic Background and Design Guidelines.....               | 135 |
| 5.1 Acoustic Background.....                                   | 135 |
| 5.1.1 Airborne.....                                            | 135 |
| 5.1.2 Structureborne.....                                      | 135 |
| 5.2 Acoustic Design Guidelines.....                            | 135 |
| 5.2.1 Magnet Room.....                                         | 135 |
| 5.2.2 Inter-Spatial Areas.....                                 | 136 |
| 5.2.2.1 Wall Construction.....                                 | 136 |
| 5.2.2.2 High Bay RF Room.....                                  | 137 |
| 5.2.2.3 Miscellaneous Plumbing, RF Windows and RF Doors.....   | 138 |
| 6 Sample Calculation AC Power Equipment Minimum Distance.....  | 139 |
| 7 Selecting Magnet Anchor Size.....                            | 141 |
| 8 Magnet Cryogenic Venting Pressure Drop Reference Tables..... | 143 |



# Chapter 1 Introduction

## 1 Preinstall Manual Introduction



### **WARNING**

PERSONNEL INJURY OR EQUIPMENT FAILURE  
FAILURE TO IMPLEMENT ALL REQUIREMENTS AND ADHERE TO ALL SPECIFICATIONS IN THIS MANUAL MAY RESULT IN PERSONAL INJURY, EQUIPMENT DAMAGE, SCAN FAILURE, OR WARRANTY VOID.

THE IMPLEMENTATION OF ALL REQUIREMENTS AND ADHERENCE TO ALL SPECIFICATIONS IN THIS MANUAL IS THE RESPONSIBILITY OF THE CUSTOMER OR THEIR ARCHITECT AND ENGINEERS. REFER ANY QUESTIONS TO THE GE HEALTHCARE PROJECT MANAGER OF INSTALLATION (PMI).



### **NOTICE**

The site must comply with all local and National codes and regulations.

### 1.1 Document Purpose

This preinstallation manual provides the necessary information to prepare a site for system installation. Specifically, this manual provides information:

1. For the site to meet system requirements
2. For the effective arrangement and interconnection of system components

### 1.2 Intended User

The primary user of this manual is the customer's installation or architectural planner who has knowledge of the following:

1. MR nomenclature, system functions, and general characteristics
2. National and local building codes
3. Customer site procedures (medical, MR, safety, etc.)
4. Any special architectural requirements (e.g., seismic codes)

### 1.3 Document Overview

This manual describes requirements and specifications for the following:

1. General System Requirements
2. Shipping and Delivery
3. Magnet Room
4. Equipment Room

- 5. Control Room
- 6. Interconnects

# Chapter 2 General System Level

## 1 Upgrade Requirements

When planning for the installation the Optima MR450w system in an existing GE Healthcare MR suite, all requirements in this manual must be met. Ensure the following:

**NOTE:** When planning installation in a non-GE Healthcare MR suite all requirements in this manual must be met.

1. The magnetic field stability tests (HSS tool) can be used to assess the vibration environment. HSS tool is currently used by local GE Field Service to evaluate the existing system performance and image quality:
  - a. Sites which currently pass system stability criteria (HSS) may proceed with the upgrade
  - b. Sites with marginal stability performance require vibration measurements to identify vibration sources
2. Remove, cover, or fill-in abandoned ducts or troughs and remove access/computer room flooring from the Equipment and Magnet rooms
3. The VibroAcoustic dampening kit must be surface mounted
4. Interconnect Cable Routing Mechanisms (ICRMs) must have a clear path in the Magnet Room (see [Chapter 6, MR System Interconnects Routing Requirements](#)). Move HVAC duct, sprinkler pipe, lights, etc. if necessary
5. RF vendor responsibilities:
  - a. The dock must be reset after the magnet (with enclosures) and patient table are aligned
  - b. The RF shield ceiling must support the cable routing mechanism. Reinforce RF shield ceiling if necessary
  - c. Dielectric breaks are required to support the ICRM
  - d. Penetration panel openings must meet the requirements in: [Chapter 4, PEN and SPW Wall Opening Requirements](#) (includes PEN and SPW panels)
  - e. RF shield attenuation must comply with: [Chapter 3, RF Shielded Room Requirements](#)
6. Install cable trays in equipment room (see [Chapter 6, MR System Interconnects Routing Requirements](#))
7. Cryogen vent must be relocated to align with the Magnet Cryogen Vent (see [Chapter 3, Magnet Room Venting Requirements](#))

## 2 System Components

The Optima MR450w system consists of the following components:

### 2.1 Magnet Room

1. 1.5T Magnet and Magnet Enclosure (MAG)
2. Rear Pedestal (PED)
3. Patient Transport Table (PT)
4. Magnet Rundown Unit (MRU). Note: An optional remote MRU may be located outside the magnet room.

### 2.2 Equipment Room

1. Main Disconnect Panel (MDP)
2. Power, Gradient, RF Cabinet (PGR)
3. Heat Exchanger Cabinet (HEC)
4. Penetration Panel Cabinet (PEN)
5. Secondary Penetration Wall (SPW)
6. Cryocooler Compressor Cabinet (CRY)
7. Magnet Monitor (MON)
8. Optional: Brainwave Lite (BW)
9. Optional: CADstream
10. Optional: MR Guided Focus Ultrasound (FUS)
11. Optional: Magnetic Resonance Elastography (MRE)

### 2.3 Control Room

1. Operator Workspace equipment (OW)
2. Pneumatic Patient Alert System (PA1)
3. Optional: Oxygen Monitor (OXY)

### 2.4 Accessories

1. Patient accessories, including phantoms, cushions, sponges, straps, and wedges
2. Gating accessories, including patient cardiac leads, peripheral gating probe, and respiratory bellows

### 3 MR Suite Minimum Room Size Requirements

Room dimensions shown in the table below lists the minimum finished room space requirements to properly and safely operate and service the MR system. The items listed below are not included in the minimum area dimensions:

1. Building code requirements (e.g., exit routes, door placement, local and national electrical codes, etc.)


**NOTE:** The customer must provide Equipment and Magnet Room evacuation routes to comply with facility emergency procedures.

2. System requirements, including cable run locations, cryogen venting, patient observation requirements, and penetration panel placements (e.g., the Equipment room and Magnet room must share a common wall to allow penetration panel installation)
3. Penetration closet (must be outside the minimum finished room dimensions)
4. Non-GEHC equipment options (such as additional AC or water cooling equipment in the Equipment room)
5. Accessory storage. Refer to the *Customer Site Storage Requirements* manual (document number 5182674) or contact the GE Healthcare Project Manager of Installation (PMI) for any additional accessory storage requirements

**Table 2-1: Minimum Finished Room Dimensions**


| Configuration                          | Equipment Room                                                                   |                            | Magnet Room<br>(See the illustrations below for specific dimensions) |                                                      |                                            | Control Room          |                                           | Total System Area<br>ft <sup>2</sup> (m <sup>2</sup> ) |
|----------------------------------------|----------------------------------------------------------------------------------|----------------------------|----------------------------------------------------------------------|------------------------------------------------------|--------------------------------------------|-----------------------|-------------------------------------------|--------------------------------------------------------|
|                                        | Area<br>ft <sup>2</sup> (m <sup>2</sup> )                                        | Ceiling Height<br>in. (mm) | W x D<br>in. (mm)                                                    | Area<br>ft <sup>2</sup> (m <sup>2</sup> )            | Finished Ceiling Height<br>in. (mm)        | W x D<br>in. (mm)     | Area<br>ft <sup>2</sup> (m <sup>2</sup> ) |                                                        |
| Minimum Room Size (2050 mm Scan Range) | 116.68 (10.84)                                                                   | 114 (2896)                 | 143.5 x 245 (3620 x 6223)                                            | Option 1: 222.49 (20.67)<br>Option 2: 218.72 (20.32) | 105 (2667)<br>Minimum Ceiling: 98.5 (2500) | 60 x 84 (1524 x 2134) | 35 (3.2)                                  | 369.85 (34.36)                                         |
| With MR Guided Focused Ultrasound      | Contact the Project Manager of Installation (PMI) for FUS site planning details. |                            |                                                                      |                                                      |                                            |                       |                                           |                                                        |
| With CADStream                         |                                                                                  |                            |                                                                      |                                                      |                                            |                       | Add: 1.25 (0.11)                          |                                                        |
| With Brainwave                         | Add: 13.25 (1.25)                                                                |                            |                                                                      |                                                      |                                            |                       |                                           |                                                        |

Illustration 2-1: Minimum Magnet Service Area

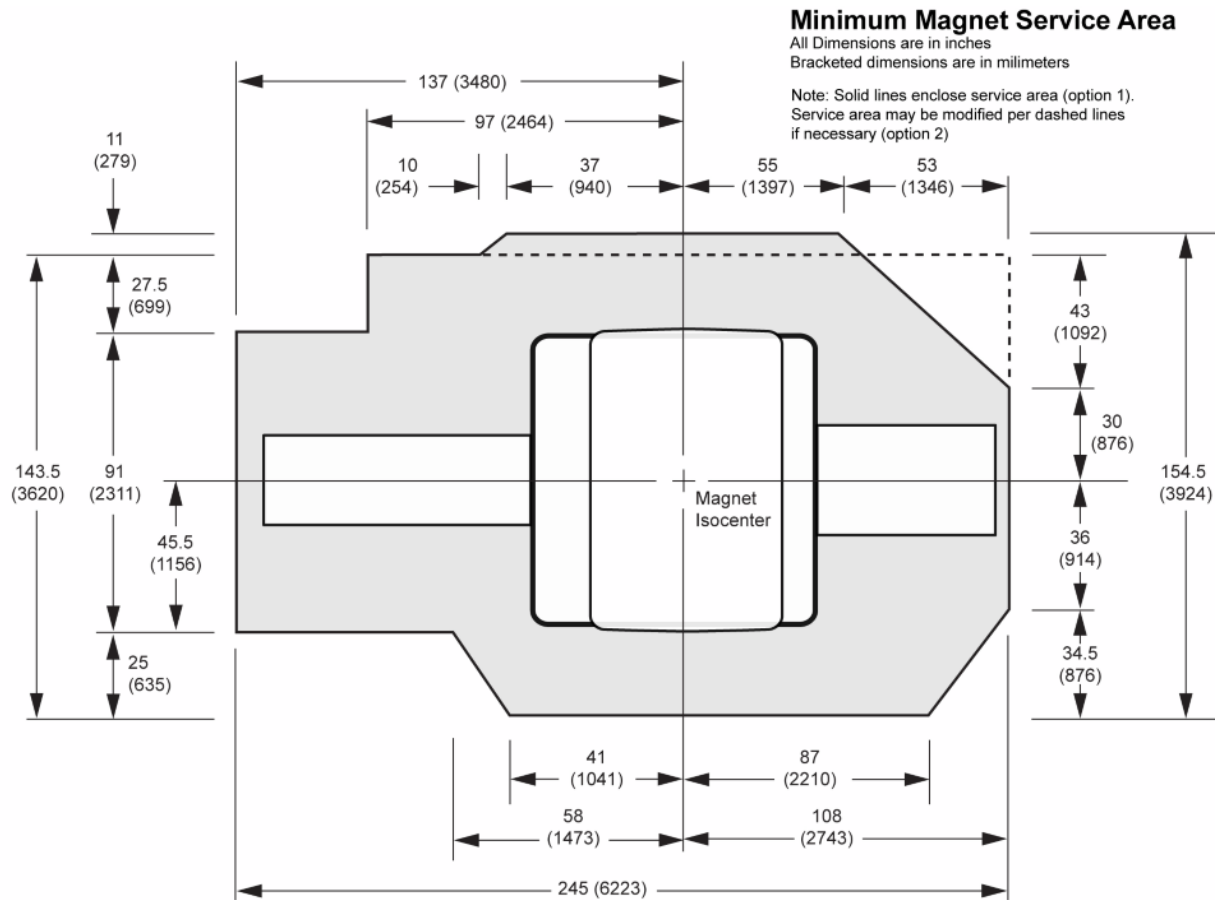**NOTICE**

The minimum ceiling height dimensions shown below are for service area only. Overall ceiling height must also take into account cryogen vent requirements [Chapter 3, Magnet Room Venting Requirements](#) and Interconnect routing requirements [Chapter 3, MR System Interconnects Routing Requirements](#).

Illustration 2-2: Minimum Magnet Ceiling Height

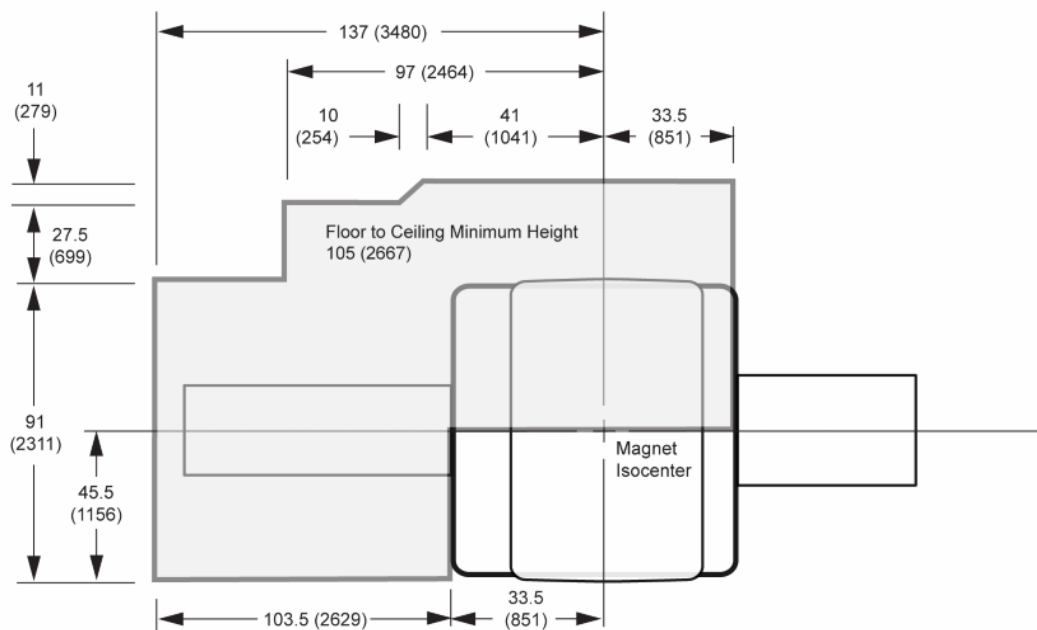**Minimum Ceiling Height**All Dimensions are in inches  
Bracketed dimensions are in millimeters

**NOTE:** Contact the GE Healthcare Project Manager of Installation (PMI) for additional information or to request designs.

## 4 IEC EMC Compliance

Per IEC 60601-1-2 Edition 2.1 Medical Electrical Equipment requires special precautions regarding Electromagnetic Compatibility (EMC) and must be installed and put into service according to the EMC information provided in the following tables. Full declaration is stored on-site in the user manual delivered with the system.

The MR system is designed and tested to the following standards:

**Table 2-2: Guidance And Manufacturer's Declaration – Electromagnetic Emissions**

| The system is intended for use in the electromagnetic environment specified below. The customer or the user of the system should assure that it is used in such an environment. |            |                                                                                                                                                                                                         |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Emissions Test                                                                                                                                                                  | Compliance | Electromagnetic Environment – Guidance                                                                                                                                                                  |
| RF Emissions CISPR 11                                                                                                                                                           | Group 2    | The system must emit electromagnetic energy in order to perform its intended function. Nearby electronic equipment may be affected                                                                      |
| RF Emissions CISPR 11                                                                                                                                                           | Class A    | The system is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes |

**Table 2-3: Guidance And Manufacturer's Declaration – Electromagnetic Immunity**

| The system is intended for use in the electromagnetic environment specified below. The customer or the user of the system should assure that it is used in such an environment. |                                                               |                                                               |
|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------|---------------------------------------------------------------|
| Immunity test                                                                                                                                                                   | IEC 60601 test level                                          | Compliance Level                                              |
| Electrostatic discharge (ESD) IEC 61000-4-2                                                                                                                                     | ±6 kV contact                                                 | ±6 kV contact                                                 |
|                                                                                                                                                                                 | ±8 kV air                                                     | ±8 kV air                                                     |
| Electrical fast transient / burst IEC 61000-4-4                                                                                                                                 | ±2 kV for power supply lines                                  | ±2 kV for power supply lines                                  |
|                                                                                                                                                                                 | ±1 kV for input/output lines                                  | ±1 kV for input/output lines                                  |
| Surge IEC 61000-4-5                                                                                                                                                             | ±1 kV line(s) to line(s)                                      | ±1 kV differential mode                                       |
|                                                                                                                                                                                 | ±2 kV line(s) to earth                                        | ±2 kV common mode                                             |
| Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11                                                                             | <5 % U <sub>T</sub> (>95 % dip in U <sub>T</sub> ) for 5 sec. | <5 % U <sub>T</sub> (>95 % dip in U <sub>T</sub> ) for 5 sec. |
| Power Frequency (50/60Hz) magnetic field IEC 61000-4-8                                                                                                                          | 3 A/m                                                         | 3 A/m                                                         |
| Conducted RF IEC 61000-4-6                                                                                                                                                      | 3 Vrms 150 kHz to 80 MHz                                      | 3 Vrms                                                        |
| Radiated RF IEC 61000-4-3                                                                                                                                                       | 3 V/m 80 MHz to 2,5 GHz                                       | 3 V/m                                                         |



## **5 MR System Seismic Requirements**

Contact the Project Manager of Installation for any questions regarding MR system seismic requirements or specifications.

1. The customer is responsible for seismic anchoring of GE components
2. Center of gravity, weight, and physical dimensions are provided for seismic calculations. Refer to the specifications or illustrations for each component

## 6 MR Suite Acoustic Specifications

### 6.1 Acoustic Specifications

The following table lists the acoustic output of GE Healthcare equipment:

**Table 2-4: Acoustic Specifications**

|                | GE Equipment Acoustic Output                                | Notes                                                                                                  |
|----------------|-------------------------------------------------------------|--------------------------------------------------------------------------------------------------------|
| Control Room   | 62 dBA                                                      |                                                                                                        |
| Equipment Room | 80 dBA                                                      | The 80 dBA level is for GE equipment only.<br>The Equipment room acoustic level must not exceed 85 dBA |
| Magnet Room    | See <a href="#">Chapter 3, Acoustic Room Specifications</a> |                                                                                                        |

Refer to [Chapter 7, Acoustic Background and Design Guidelines](#)

**NOTE:** All GE equipment acoustic output values are for base equipment configuration in each room.

### 6.2 Structureborne Vibration Control Specifications

Structureborne acoustic issues tend to occur at MR installations above the ground floor of the facility. Two options to mitigate structureborne acoustic transmission are:

1. GE Healthcare provides a VibroAcoustic Dampening kit (which must be surface mounted). Contact the GE Healthcare Project Manager of Installation for information
2. The customer may design and implement a custom solution in addition to the VibroAcoustic Dampening kit

**NOTE:** The amount of vibration attenuation provided by the VibroAcoustic Dampening kit will be site dependent.

## 7 MR Suite Magnetic Field Specifications

### 7.1 Magnetic Fringe Field

The following illustrations show the static magnet isogauss plot lines for a 1.5T magnet. This information must be used to evaluate potential site interaction with internal and external magnetic fields, magnetic materials on the site, and to locate personnel and equipment within the site.

The 5 gauss line can expand to 14.27 ft (4.35 m) axially and 9.35 ft (2.85 m) radially for up to 1 second in the rare event of a quench

These isogauss plots show an idealized magnetic field relative to magnet isocenter. The actual field strength can be affected by any of the following:

- Magnetic shielding
- Earth's magnetic field
- Other magnetic fields
- Stationary or moving metal

Illustration 2-3: Magnetic Fringe Field Side View

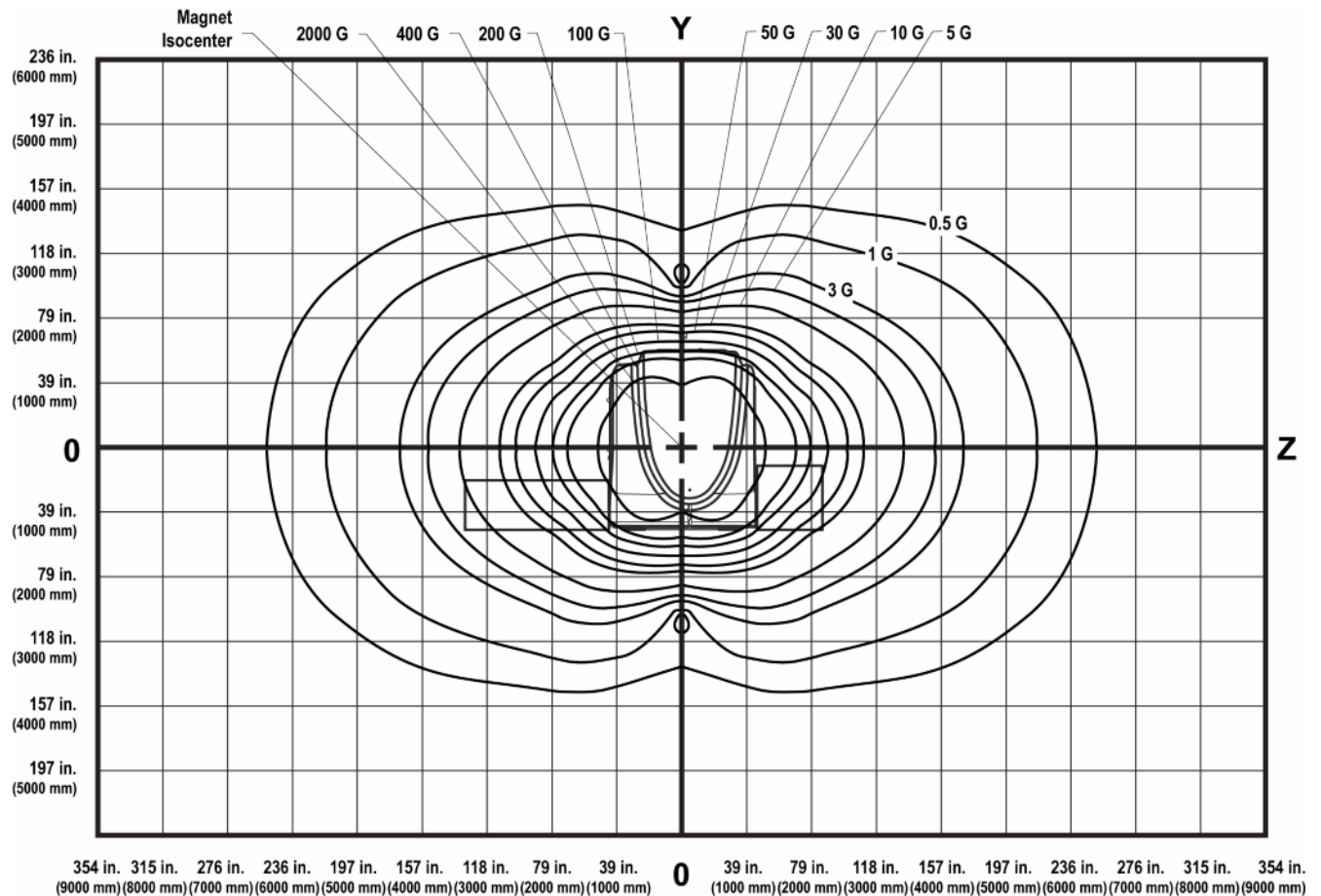


Illustration 2-4: Magnetic Fringe Field Top View


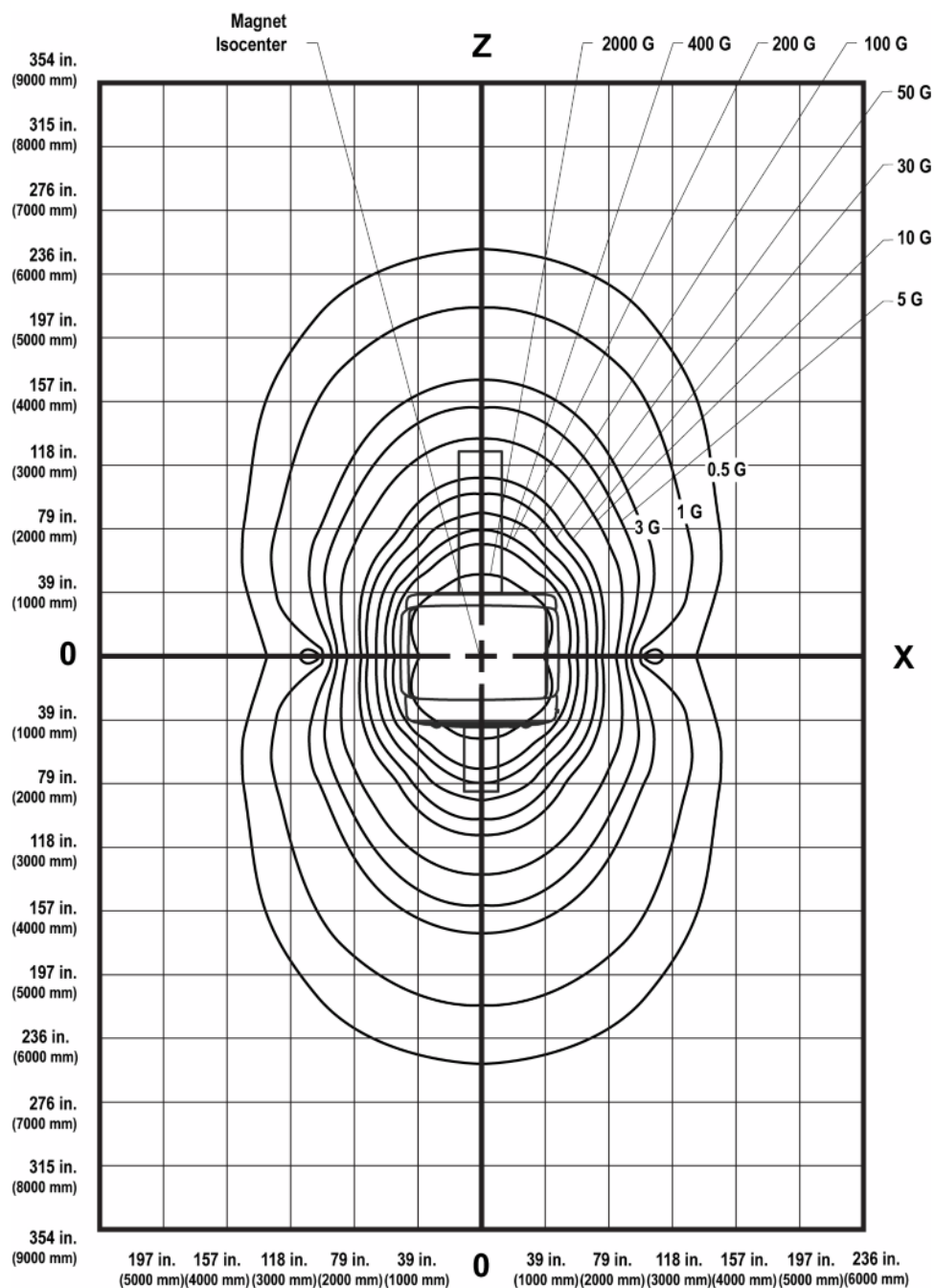


Illustration 2-5: Magnetic Fringe Field Front View

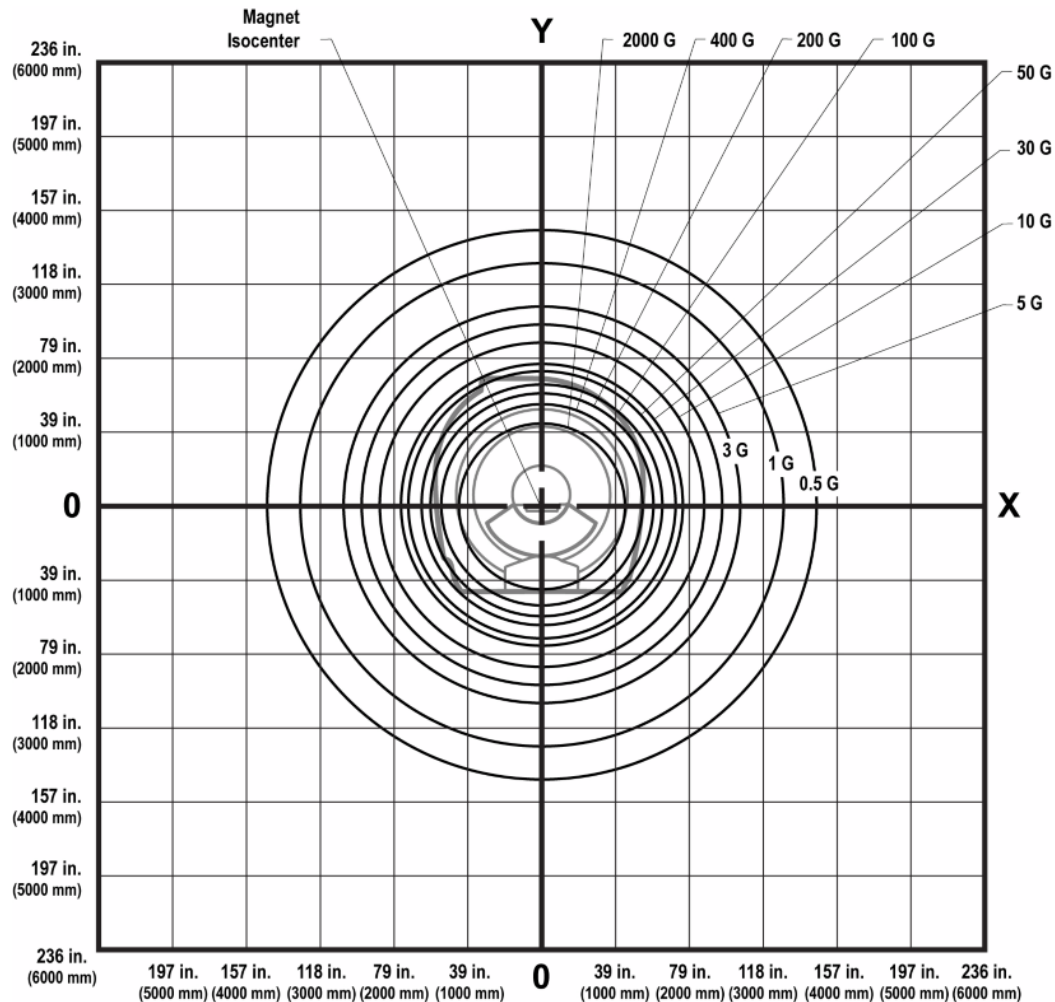

## 7.2 Interference from Changing Magnetic Fields

Metal objects moving within the magnet sensitivity lines can produce a field disturbance during clinical imaging. If the metal object is moving it will produce a fluctuating dipole type of field which cause image artifacts. As an example, a car driven inside the moving metal line will act as a dipole and produce a time varying field which change the magnet's main field during the imaging time. The same vehicle may park within the moving metal line and remain parked during clinical scanning without impact to the main field. See [Illustration 2-6](#) and [Illustration 2-7](#).

Illustration 2-6: Magnet Moving Metal Sensitivity Line Plot (Top View)

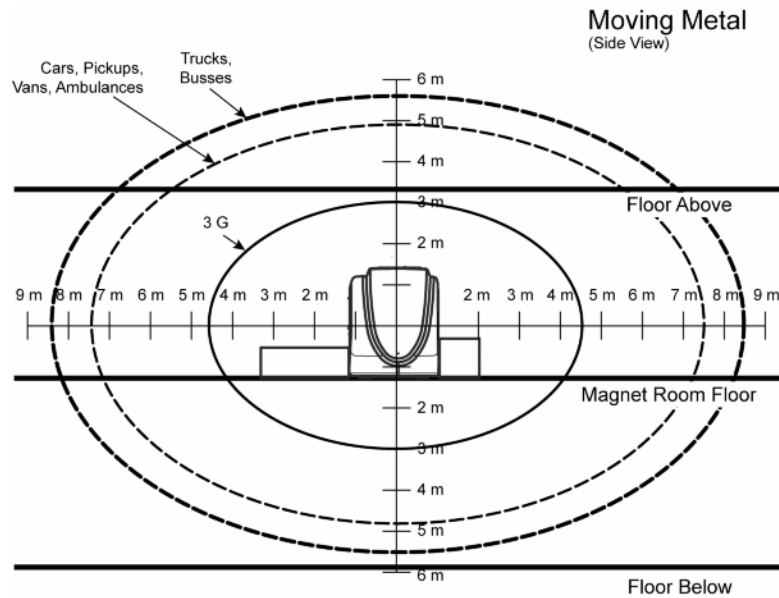


Illustration 2-7: Magnet Moving Metal Sensitivity Line Plot (Side View)

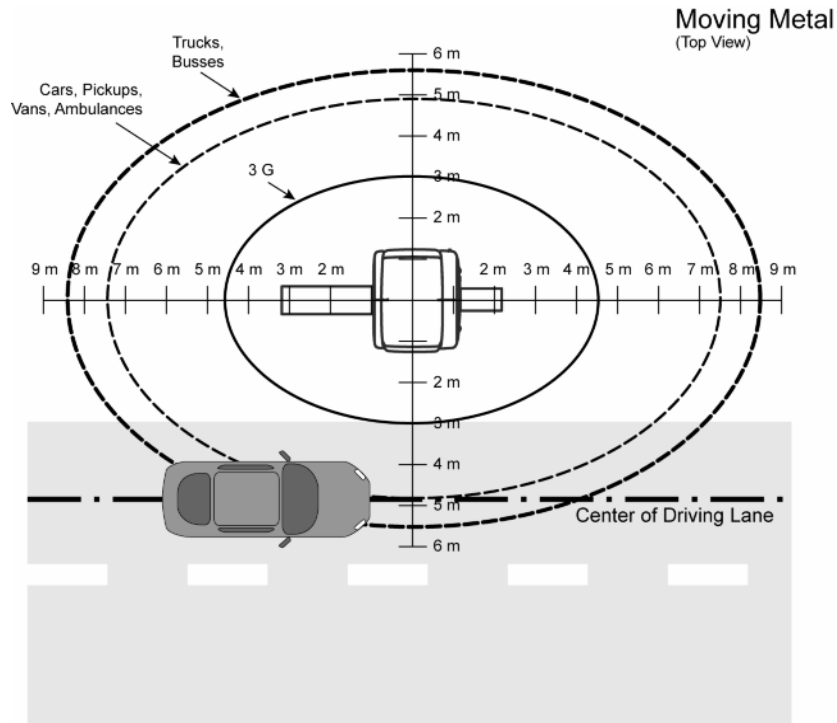
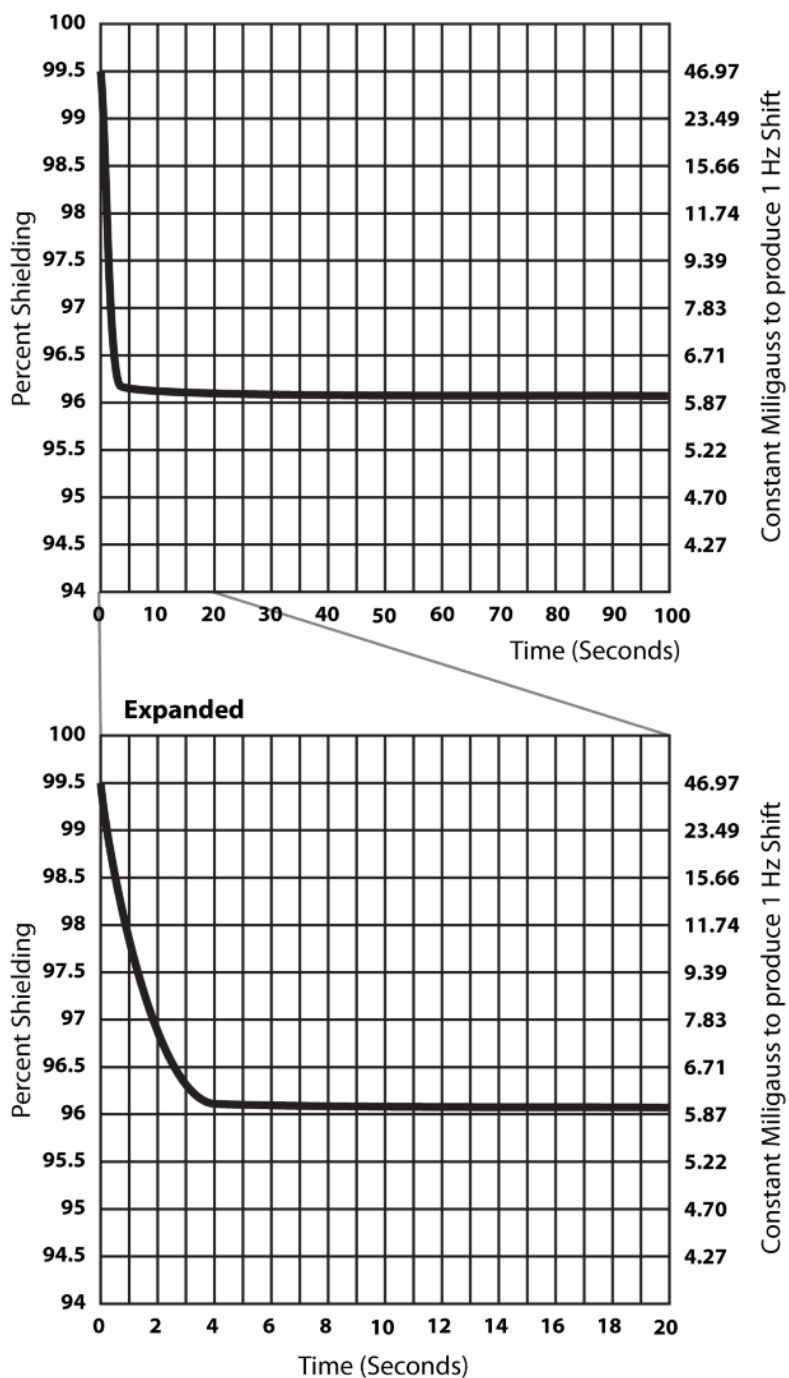


Table 2-5: Magnet Moving Metal Requirements

| Metal Objects Category | Definition Of Distance Location        | Magnet Minimum Distance Radial X Axial ft (m) See Note 1 |
|------------------------|----------------------------------------|----------------------------------------------------------|
| Objects 100 - 400 lbs  | Distance from isocenter radial x axial | 3 Gauss line                                             |

| Metal Objects Category                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         | Definition Of Distance Location                                                                                                                                       | Magnet Minimum Distance Radial X Axial ft (m) See Note 1   |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------|
| Cars, Minivans, Vans, Pickup Trucks, Ambulances                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                | Distance from isocenter measured to center of driving or parking lane radial x axial                                                                                  | 15.5 x 24.6 (4.72 x 7.5)                                   |
| Bus, Trucks (Utility, Dump, Semi)                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                              | Distance from isocenter measured to center of driving or parking lane radial x axial                                                                                  | 18.1 x 28.75 (5.52 x 8.76)                                 |
| Objects > 400 lbs, Elevators, Trains, Subways                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  | Place a directional probe (e.g. flux gate sensor) at isocenter of proposed magnet location aligned along the Z-axis. Measure peak-to-peak magnetic field change (DC). | Refer to <a href="#">Illustration 2-8</a> and notes below: |
| <p>Notes:</p> <ol style="list-style-type: none"> <li>Radial distances are magnet X and Y axis. Axial distances are magnet Z axis.</li> <li>EXAMPLE: For Moving Metal Requirements of objects &gt; 400 lbs category you can use the time history of the occurrence to determine what milligauss level to use. <ol style="list-style-type: none"> <li>If the site has elevators/counter weights near the magnet and the elevator can stop on the floors for longer than 20 seconds (which is usually the case), peak-to-peak milligauss reading must be less than 5.87.</li> <li>If the site has a subway nearby and the field disturbance is less than 2 seconds, the peak-to-peak milligauss reading must be less than 5.87.</li> <li>Use 5.87 milligauss peak-to-peak.</li> </ol> </li> </ol> |                                                                                                                                                                       |                                                            |

Illustration 2-8: Actual Axial Shielding Performance

Actual Axial  
Shielding Performance

### 7.3 Electrical Current


1. Electrical current in high voltage power lines, transformers, motors, or generators near the magnet may affect magnetic field homogeneity



2. Magnetic field interference at 50 or 60 Hz must not exceed 40 milligauss RMS respectively at the magnet location (refer to [Illustration 2-9](#))
3. The following equation can be used as a general guide in determining allowable current in feeder lines at a given distance from the magnet isocenter:
  - a. For 1.5T Magnet:  $I = (20X^2)/S$
  - b. I = Maximum allowable RMS single phase current (in amps) or maximum allowable RMS line current (in amps) in three phase feeder lines
  - c. S = Separation (in meters) between single phase conductors or greatest separation between three phase conductors
  - d. X = Minimum distance (in meters) from the feeder lines to isocenter of the magnet

**Illustration 2-9: Magnet Allowable Milligauss vs. Line Frequency for AC Equipment**

AC Field EMI Limits

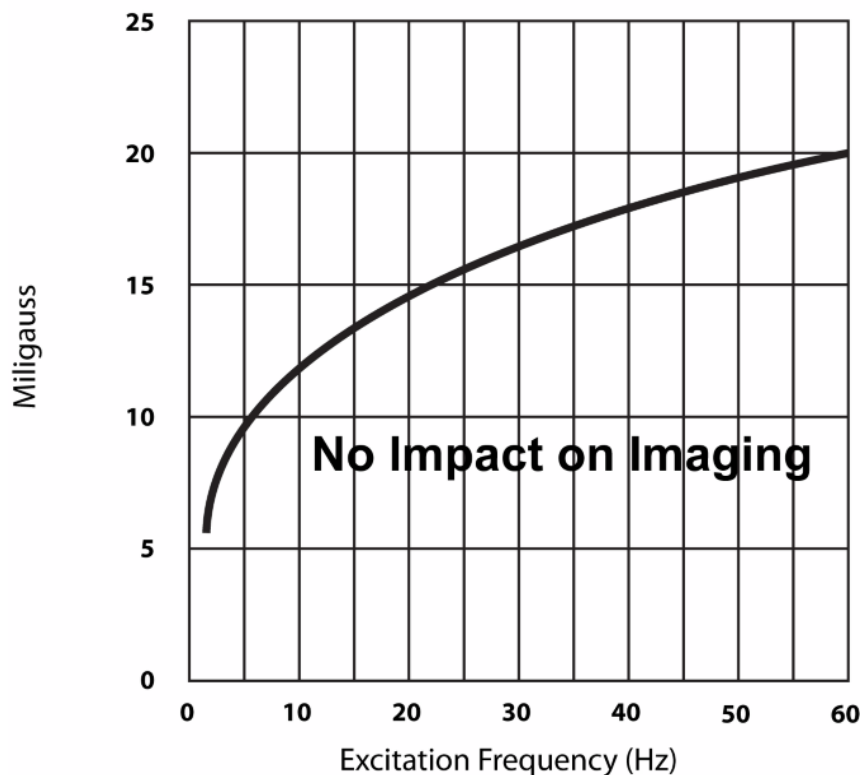

Refer to [Chapter 7, Sample Calculation AC Power Equipment Minimum Distance](#) for additional examples.

## 7.4 Non-MR System Equipment Sensitivity to Magnetic Fields

This section lists equipment known to be sensitive to high magnetic fields. Recommended limits given are based on general MR site planning guidelines. Actual susceptibility of specific devices

may vary significantly depending on electrical design, orientation of the device relative to the magnetic field, and the degree of interference considered unacceptable.

The table is provided as a guide only. Actual Gauss limits are equipment specific (refer to the OEM manuals). Site plans must include consideration for magnetic field interaction with customer equipment.

**NOTE:** Gauss limits for MR system equipment is provided with the specifications for each component.

**Table 2-6: Magnetic Proximity Limits**


| Gauss (mT) Limit   | Equipment                            |                                           |
|--------------------|--------------------------------------|-------------------------------------------|
| 0.5 gauss (0.05mT) | Nuclear camera                       |                                           |
| 1 gauss (0.1mT)    | Positron Emission Tomography scanner | Video display (tube)                      |
|                    | Linear Accelerator                   | CT scanner                                |
|                    | Cyclotrons                           | Ultrasound                                |
|                    | Accurate measuring scale             | Lithotripter                              |
|                    | Image intensifiers                   | Electron microscope                       |
|                    | Bone Densitometers                   |                                           |
| 3 gauss (0.3mT)    | Power transformers                   | Main electrical distribution transformers |
| 5 gauss (0.5mT)    | Cardiac pacemakers                   | Biostimulation devices                    |
|                    | Neurostimulators                     |                                           |
| 10 gauss (1mT)     | Magnetic computer media              | Telephone switching stations              |
|                    | Hard copy imagers                    | Water cooling equipment                   |
|                    | Line printers                        | HVAC equipment                            |
|                    | Video Cassette Recorder (VCR)        | Major mechanical equipment room           |
|                    | Film processor                       | Credit cards, watches, and clocks         |
|                    | X-ray tubes                          |                                           |
|                    | Large steel equipment, including:    |                                           |
|                    | Emergency generators                 | Air conditioning equipment                |
|                    | Commercial laundry equipment         | Fuel storage tanks                        |
|                    | Food preparation area                | Motors greater than 5 horsepower          |
|                    |                                      |                                           |
| 50 gauss (5mT)     | Metal detector for screening         | Telephones                                |
|                    | LCD panels                           |                                           |
| No Limit           | Digital Detectors                    |                                           |

## 8 Multiple MR System Requirements

### 8.1 Multiple Magnets

When installing multiple magnets, the 3 gauss lines must not intersect. Contact the GE Healthcare Project Manager of Installation (PMI) for any questions regarding magnetic field interaction.

Illustration 2-10: Two Magnet Installation

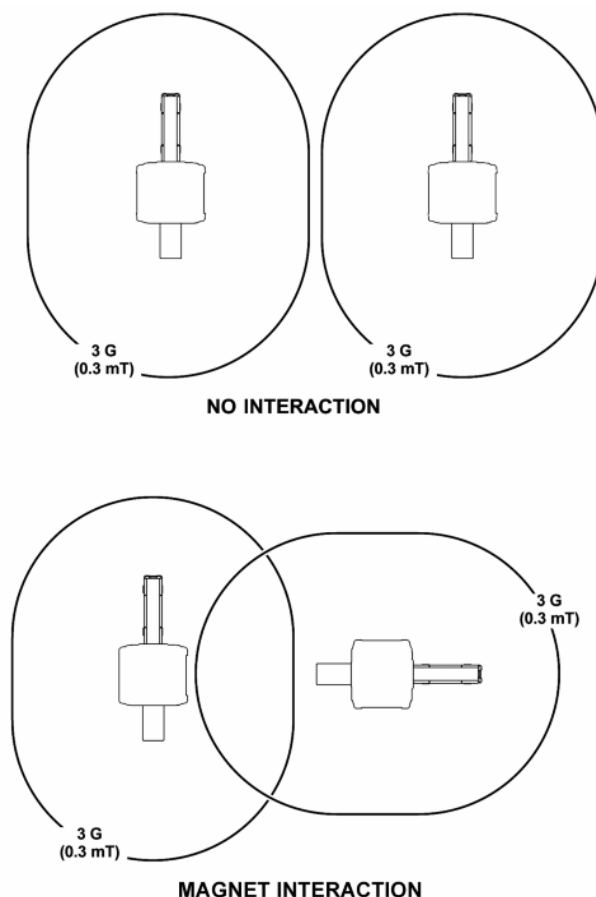

### 8.2 Shared Equipment Rooms

When installing multiple MR systems in a shared equipment room, the following conditions must be met:

1. RF cabinets must be separated by at least 79 in. (2000 mm). Refer to [Illustration 2-11](#)
2. PEN Panels must be separated by at least 118 in (3000 mm)
3. Cables from different MR systems must not share interconnect cable routing mechanisms



## Illustration 2-11: Equipment Room Cabinet Separations

## Shared Equipment Room

## Minimum Separation Distances

All Dimensions are in inches

Bracketed dimensions are in millimeters

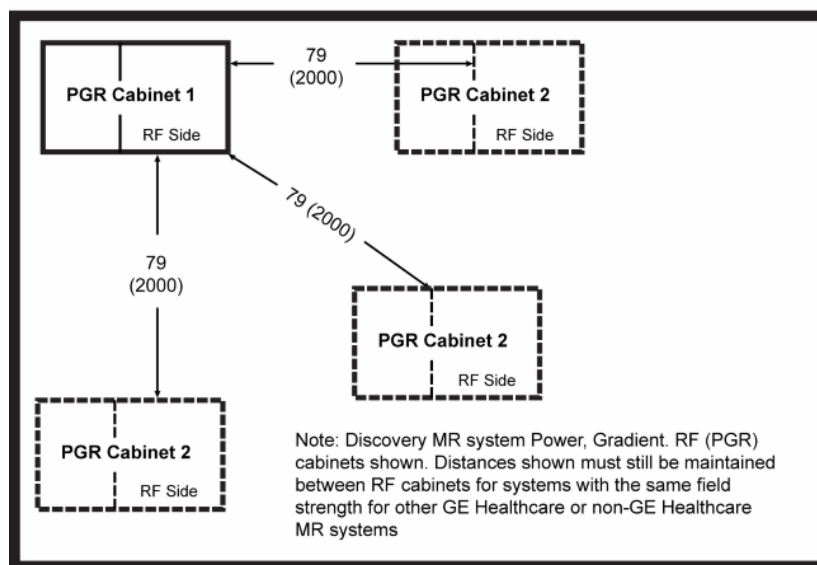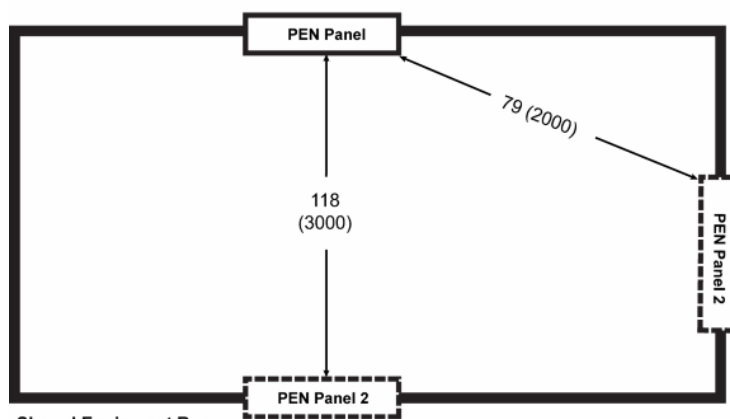Shared Equipment Room  
(Top View)Shared Equipment Room  
(Top View)

## 9 MR Suite Temperature and Humidity



### CAUTION

#### Equipment Failure

Failure to maintain the required temperature or humidity at all times (i.e., both working and non-working hours) may result in equipment failure, scanning failure, or warranty void.

Ensure the HVAC system has the correct capacity for the room size, equipment heat output, and environmental conditions to maintain proper temperature and humidity.

This section provides temperature and humidity requirements for the MR suite. Specific construction requirements for each room can be found in the following sections:

- Magnet Room
- Equipment Room
- Control Room

### 9.1 Temperature and Humidity Requirements

Table 2-7: Room Temperature and Humidity Requirements

| Room                                                                                                                           | Temperature                  |                                   | Humidity  |                            |
|--------------------------------------------------------------------------------------------------------------------------------|------------------------------|-----------------------------------|-----------|----------------------------|
|                                                                                                                                | Range °F (°C)                | Change °F/Hr (°C/Hr) <sup>1</sup> | Range %RH | Change %RH/Hr <sup>2</sup> |
| Equipment Room (at Inlet to Equipment)                                                                                         | 59-89.6 (15-32) <sup>3</sup> | 5 (3)                             | 30-70     | 5                          |
| Magnet Room                                                                                                                    | 59-69.8 (15-21)              | 5 (3)                             | 30-60     | 5                          |
| Operator Room                                                                                                                  | 59-89.6 (15-32)              | 5 (3)                             | 30-70     | 5                          |
| Notes:                                                                                                                         |                              |                                   |           |                            |
| 1. Operating temperature gradient limits shall be between -5° F/Hr (-3° C/Hr) and 5° F/Hr (3° C/Hr), when averaged over 1 hour |                              |                                   |           |                            |
| 2. Operating humidity gradient limits shall be between -5% RH/hour and 5% RH/hour, when averaged over 1 hour                   |                              |                                   |           |                            |
| 3. Maximum ambient temperature is derated by 1 degree C per 300 m above 2000 m                                                 |                              |                                   |           |                            |

1. The customer is responsible for HVAC system design, purchase, and installation
2. The temperature requirements must not be exceeded at any point during the day (both working or non-working hours)
3. A separate thermostat must be provided for the Magnet room

### 9.2 Equipment Heat Output Specifications

This section details the heat output for specific components. These heat outputs define the minimum, maximum and an assumed average condition. Actual heat output and room temperature may vary due to environmental factors, room insulation, actual usage, and any non-GE Healthcare equipment used in the MR suite. Also, due to large variations in heat loads, the HVAC system may require unloaders, hot gas bypass, and reheat to maintain humidity levels.

**NOTE:** Air cooling load for components in the Equipment room is averaged over a 12 hour working day is approximately 1/2 of the maximum value.

**Table 2-8: System Heat Output for Air Cooling**

| Component                                     | Magnet Room   |               |             | Equipment Room |                |               | Control Room  |         |      |
|-----------------------------------------------|---------------|---------------|-------------|----------------|----------------|---------------|---------------|---------|------|
|                                               | BTU/hr (W)    |               |             | BTU/hr (W)     |                |               | BTU/hr (W)    |         |      |
|                                               | Maximum       | Average       | Idle        | Maximum        | Average        | Idle          | Maximum       | Average | Idle |
| Magnet (MAG) and Patient Transport Table (PT) | 8,189 (2,400) | 4,095 (1,200) | 1,915 (561) |                |                |               |               |         |      |
| Blower Box                                    | 1,535 (450)   | 1,535 (450)   | 1,535 (450) |                |                |               |               |         |      |
| Penetration Panel Cabinet (PEN)               | 1,024 (300)   | 512 (150)     | 0           | 10,697 (3,135) | 5,349 (1,568)  | 5,349 (1,568) |               |         |      |
| Secondary Penetration Wall (SPW)              | 0             |               |             | 0              |                |               |               |         |      |
| Main Disconnect Panel (MDP)                   |               |               |             | 901 (264)      | 450 (132)      | 450 (132)     |               |         |      |
| Power, Gradient, RF Cabinet (PGR)             |               |               |             | 20,940 (6,137) | 10,470 (3,068) | 8,530 (2500)  |               |         |      |
| Cryocooler Compressor (CRY)                   |               |               |             | 1,706 (500)    | 1,706 (500)    | 1,706 (500)   |               |         |      |
| Heat Exchanger Cabinet (HEC)                  |               |               |             | 3,412 (1,000)  | 1,706 (500)    | 1,706 (500)   |               |         |      |
| Magnet Monitor (MON)                          |               |               |             | 819 (240)      | 819 (240)      | 819 (240)     |               |         |      |
| Operator Workspace equipment (OW)             |               |               |             |                |                |               | 4,947 (1,450) |         |      |

**Table 2-9: System Options Heat Output for Air Cooling**

| Component                      | Magnet Room |         |      | Equipment Room                            |         |      | Control Room |         |      |
|--------------------------------|-------------|---------|------|-------------------------------------------|---------|------|--------------|---------|------|
|                                | BTU/hr (W)  |         |      | BTU/hr (W)                                |         |      | BTU/hr (W)   |         |      |
|                                | Maximum     | Average | Idle | Maximum                                   | Average | Idle | Maximum      | Average | Idle |
| BrainWave HW Lite Cabinet (BW) |             |         |      | 2337 (685) base (2781 (815) with options) |         |      |              |         |      |
| CADstream                      |             |         |      | Contact PMI                               |         |      |              |         |      |
| MR Elastography                |             |         |      | 480                                       |         |      |              |         |      |

## 10 Facility Coolant Requirements



### NOTICE

**Equipment Failure.** A continuous supply of facility liquid coolant to the Heat Exchanger Cabinet (HEC) is required at all times for proper system operation. Failure to provide liquid coolant with the requirements listed in this section may cause equipment failure.

### 10.1 Heat Exchanger Cabinet (HEC) Coolant Requirements

1. The facility must provide liquid coolant to the Heat Exchanger Cabinet (HEC)
2. The facility must provide pipe/hose, filter, and connectors to the HEC
3. The facility must provide an uninterrupted supply of liquid coolant to the system at magnet delivery
4. The vertical distance between the coolant connection points of the HEC and the Gradient Coil must be less than 196.8 in. (5 meters)
5. The customer must provide and install an in-line flow meter on either the supply or return facility coolant hose. The flow meter must be capable of displaying volumetric flow in either gallons per minutes or liters per minute.


**Table 2-10: Facility Liquid Coolant Requirements**

| Parameter                                                     | Requirements                                                                                                            |
|---------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------|
| Availability                                                  | Continuous                                                                                                              |
| Antifreeze                                                    | 0-40% propylene glycol                                                                                                  |
| Minimum Flow                                                  | 30 gpm (114 L/min)                                                                                                      |
| Maximum Flow                                                  | 35 gpm (132 L/min)                                                                                                      |
| Maximum Pressure Drop in HEC at Minimum Flow                  | 30 psi (2.4 bar) with 40% propylene glycol-water; 1021 kg/m <sup>3</sup> density                                        |
| Maximum Pressure Drop in HEC at Maximum Flow                  | 50 psi (3.4 bar) with 40% propylene glycol-water; 1021 kg/m <sup>3</sup> density                                        |
| Temperature rise at Minimum Flow                              | 12.2°F (6.8°C) with 40% propylene glycol-water; 3730 J/(kg K) specific heat; 1021 kg/m <sup>3</sup> density; 49 kW heat |
| Temperature rise at Maximum Flow                              | 10.4°F (5.8°C) with 40% propylene glycol-water; 3730 J/(kg K) specific heat; 1021 kg/m <sup>3</sup> density; 49 kW heat |
| Inlet Pressure after the filter                               | 57 to 87 psi (3.9 to 6 bar) measured at the inlet to the HEC                                                            |
| Chiller Size                                                  | Minimum 49 KW                                                                                                           |
| Condensation Protection                                       | Condensation must be managed to prevent equipment damage or safety hazards.                                             |
| Minimum Continuous Heat Load                                  | 7.5 KW                                                                                                                  |
| Inlet Temperature                                             | 44.6 to 50°F (7 to 10°C) measured at the inlet to the HEC                                                               |
| Customer supplied feeder hose (from main water supply to HEC) | 1.5 inch (38.1 mm) minimum hose inside diameter                                                                         |
| Hose connections to the HEC                                   | 1.5 inch (38.1 mm) male NPT                                                                                             |

Table 2-11: Water Quality

|                         |                                             |
|-------------------------|---------------------------------------------|
| pH Value                | 6.5 to 8.2 at 77 °F (25 °C)                 |
| Electrical Conductivity | < 0.8 mmho/cm                               |
| Chloride Ion            | < 200 ppm                                   |
| Sulfate Ion             | < 200 ppm                                   |
| M-Alkalinity            | < 100 ppm                                   |
| Total Hardness          | < 200 ppm                                   |
| Calcium Hardness        | < 150 ppm                                   |
| Ionic Silica            | < 50 ppm                                    |
| Iron                    | < 1.0 ppm                                   |
| Copper                  | < 0.3 ppm                                   |
| Sulfide Ion             | None, not detectable                        |
| Ammonium Ion            | < 1.0 ppm                                   |
| Residual Chlorine       | < 0.3 ppm                                   |
| Free Carbon Dioxide     | < 4.0 ppm                                   |
| Stability Index         | 6.0 to 7.0                                  |
| Suspended Matter        | < 10 ppm                                    |
| Particle Size           | < 100 micron (with field changeable filter) |

Illustration 2-12: MR System Water Cooling Block Diagram

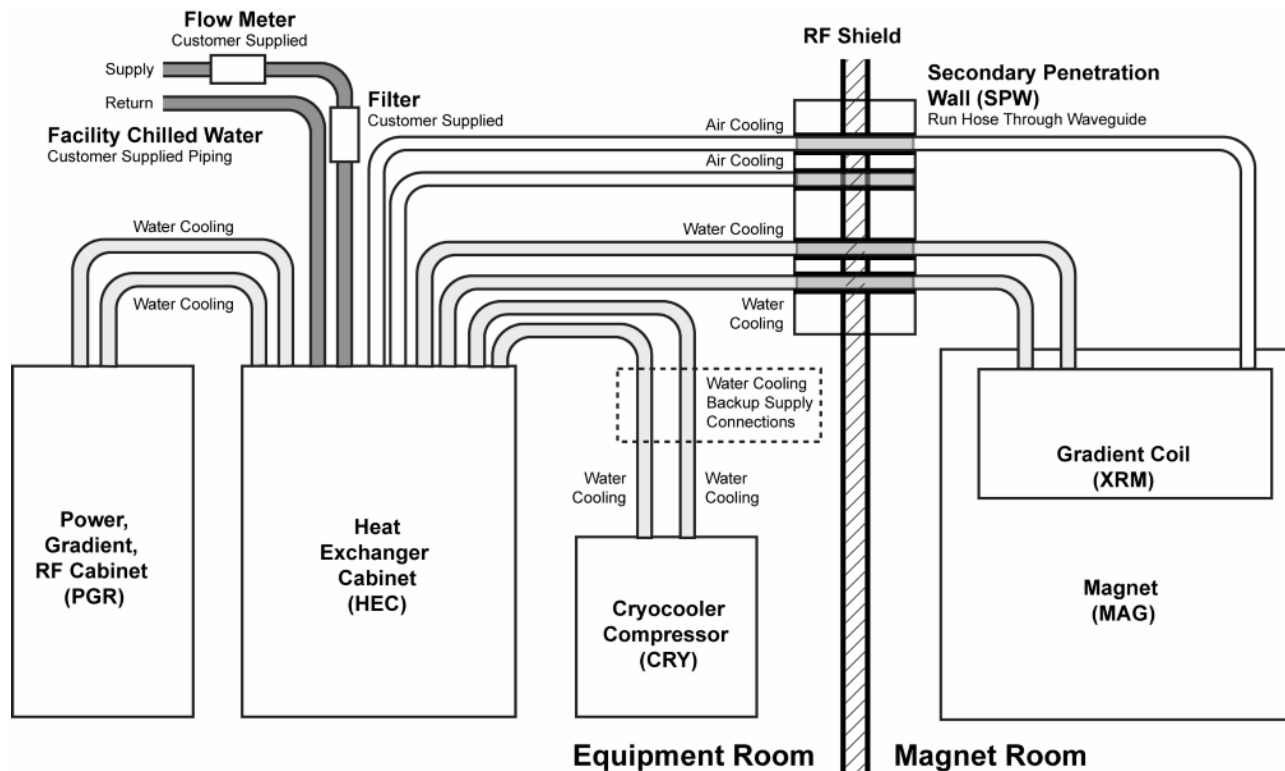

## 10.2 Cryocooler Compressor (CRY) Emergency Facility Coolant Requirements



The facility may provide an optional backup coolant supply to the Cryocooler Compressor.



1. The facility is responsible for coolant, pipe/hose, filters, and connectors to supply the coolant to the CRY
2. The emergency coolant supply must isolate the Cryocooler Compressor and not back-feed the HEC
3. Coolant must meet all other HEC coolant specifications listed in [Table 2-10](#) and [Table 2-11](#)
4. The supplied water cooling hoses between the HEC and CRY are Parker Push-lok 801-8, 0.5 in. I.D.
5. The charts below shows the coolant flow rate and temperature requirements for the Cryocooler Compressor:



Illustration 2-13: Cryocooler Water Cooling requirements for Emergency Water Supply

**For Water**

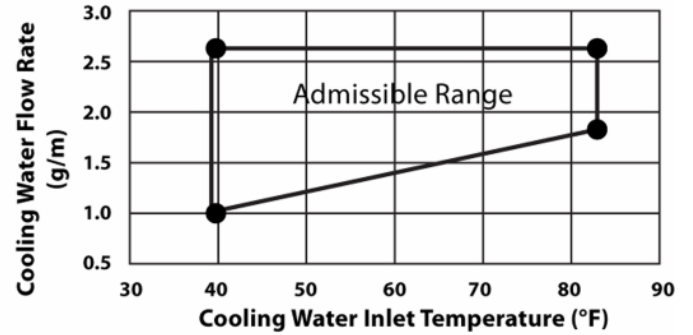
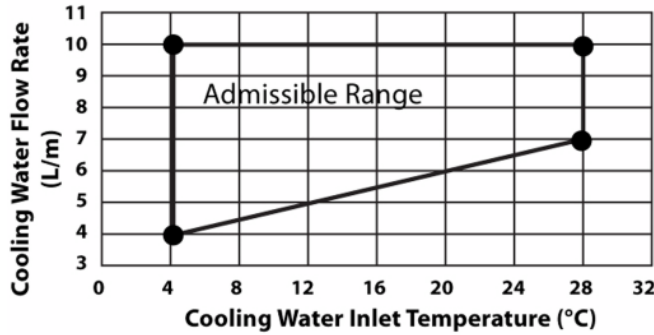


**For Antifreeze**

(50/50% mixture of water and ethylene glycol)

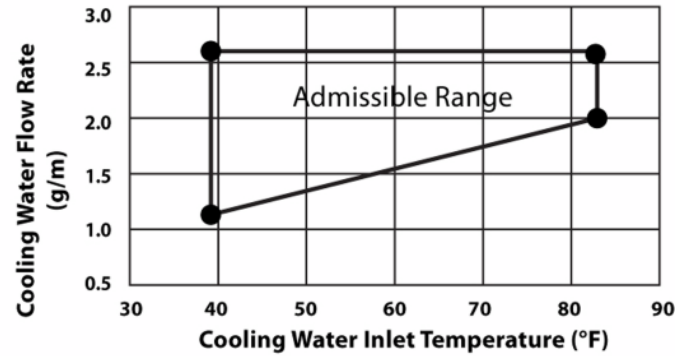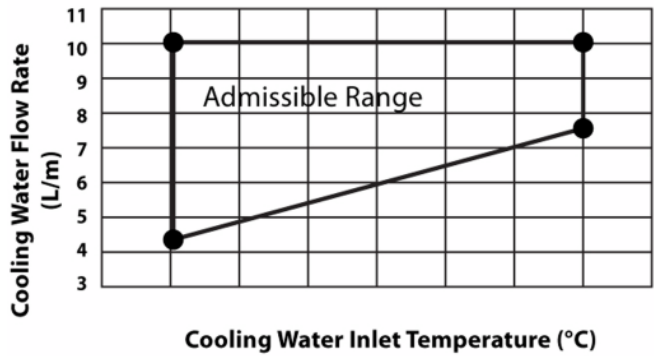

**Pressure Drop and Temperature Rise**

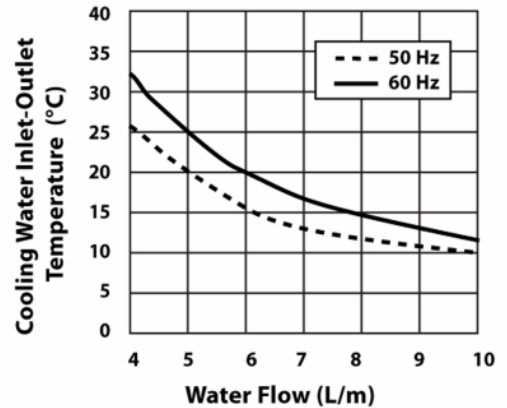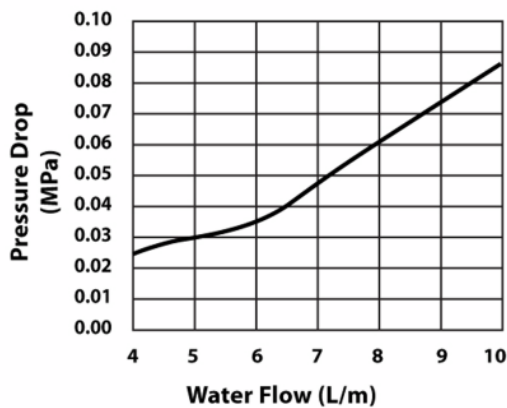

## 11 MR Suite Electrical Requirements

1. An MDP is provided with the system and is delivered prior to system delivery
2. The facility must provide system power and cabling to the Main Disconnect Panel (MDP)
3. All associated transformers and cables must be correctly sized for system power requirements
4. Facility power is required at magnet delivery and must be available continuously thereafter


**Table 2-12: Facility Power Requirements**

| Component                   | Parameter                              | Requirements                                                                                                                                                                                                                                                                                                                                                                                                                                                                                |                                                                    |
|-----------------------------|----------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------|
| Main Disconnect Panel (MDP) | Voltage / Frequency                    | 480 VAC                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 60 $\pm$ 3 Hz                                                      |
|                             |                                        | 415 VAC                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 50/60 $\pm$ 3 Hz                                                   |
|                             |                                        | 400 VAC                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 50/60 $\pm$ 3 Hz                                                   |
|                             |                                        | 380 VAC                                                                                                                                                                                                                                                                                                                                                                                                                                                                                     | 50/60 $\pm$ 3 Hz                                                   |
|                             | Daily Voltage Variation                | $\pm$ 10% from nominal under worst-case line and load regulation                                                                                                                                                                                                                                                                                                                                                                                                                            |                                                                    |
|                             | Voltage Transients                     | Maximum allowable transient voltage above or below nominal waveshape not to exceed 200 V at a maximum duration of 1 cycle and frequency of 10 times per hour.                                                                                                                                                                                                                                                                                                                               |                                                                    |
|                             | Phase                                  | Input power to the MDP may use one of the following configurations: <ul style="list-style-type: none"> <li>• A 3 phase solidly grounded WYE with Neutral and Ground (5-wire system). Neutral must be terminated prior to or inside the Main Disconnect Panel and not brought to the Power, Gradient, Reconstruction cabinet (PGR) or the Heat Exchanger Cabinet (HEC).</li> <li>• A 3 phase DELTA with Ground (4-wire).</li> </ul>                                                          |                                                                    |
|                             | Phase Balance                          | Difference between the highest phase line-to-line voltage and the lowest phase line-to-line voltage must not exceed 2%                                                                                                                                                                                                                                                                                                                                                                      |                                                                    |
|                             | Regulation                             | 2% maximum at system maximum power demand (averaged over 5 seconds) from source to PDU input (i.e. includes all feeders and transformer to utility)                                                                                                                                                                                                                                                                                                                                         |                                                                    |
|                             | Facility Zero Voltage Reference Ground | <ul style="list-style-type: none"> <li>• The facility ground for the MR system must originate at the system power source (i.e., transformer or first access point of power into the facility) and be continuous to the MR system Main Disconnect Panel (MDP) in the room.</li> <li>• Main facility ground conductor to Main Disconnect Panel (MDP) must be appropriately sized insulated copper wire.</li> <li>• Ground impedance to earth at power source to be 2 ohms or less.</li> </ul> |                                                                    |
|                             | Power Demand                           | Maximum Momentary                                                                                                                                                                                                                                                                                                                                                                                                                                                                           | 123 kVA for 5 sec consisting of 103 kVA for PDU and 20 kVA for HEC |
|                             |                                        | Maximum Continuous (while scanning)                                                                                                                                                                                                                                                                                                                                                                                                                                                         | 99 kVA consisting of 79 kVA for PDU and 20 kVA for HEC             |
|                             |                                        | Standby (no scan)                                                                                                                                                                                                                                                                                                                                                                                                                                                                           | No more than 17 kVA                                                |
|                             | Power Availability                     | Continuous. Facility power is required at all times for operation of the Heat Exchanger Cabinet (HEC) and Cryocooler (CRY).                                                                                                                                                                                                                                                                                                                                                                 |                                                                    |
|                             | MDP Circuit Breaker                    | CB#1: 200 Amps, CB#2: 150 Amps, CB#3: 50 Amps                                                                                                                                                                                                                                                                                                                                                                                                                                               |                                                                    |
|                             | MDP Short Circuit Current Rating       | 25,000 Amps                                                                                                                                                                                                                                                                                                                                                                                                                                                                                 |                                                                    |
|                             | Facility Supplied Wiring               | Facility supplied wiring is required for facility input power to MDP, from MDP to PGR, and from MDP to HEC.                                                                                                                                                                                                                                                                                                                                                                                 |                                                                    |
|                             |                                        | MDP phase input CB#1                                                                                                                                                                                                                                                                                                                                                                                                                                                                        | Accommodates 3/0 AWG maximum wire size                             |
|                             |                                        | MDP Input Ground/PE terminal                                                                                                                                                                                                                                                                                                                                                                                                                                                                | 8AWG to 350KCMIL                                                   |
|                             |                                        | MDP Input Neutral Terminal                                                                                                                                                                                                                                                                                                                                                                                                                                                                  | Accommodates 3/0 AWG maximum wire size                             |

| Component                         | Parameter                                                                      | Requirements             |                                                                                                                                                                             |
|-----------------------------------|--------------------------------------------------------------------------------|--------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|                                   |                                                                                | PGR Input phase & ground | 1/0 AWG wire size                                                                                                                                                           |
|                                   |                                                                                | HEC Input phase & ground | 8 AWG wire size (or equivalent). Provided strain relief accepts non-metallic jacketed cord with outside diameter range: 0.875 in. (22.0 mm) min. to 1.26 in. (32.0 mm) max. |
| Emergency Off Buttons             | Two Remote flush wall mounted EMERGENCY OFF buttons are supplied with the MDP. |                          |                                                                                                                                                                             |
| Service receptacle in Magnet Room | Voltage / Frequency                                                            | 100-120 VAC 50/60 Hz     | Receptacle required for small power tools. Local voltage and portable transformers for voltage values.                                                                      |
|                                   | Phase                                                                          | 1                        |                                                                                                                                                                             |
|                                   | Maximum Amps                                                                   | 2.0                      |                                                                                                                                                                             |

Illustration 2-14: MR System Main Disconnect Panel (MDP) Set-Up

- Note:**
- Runs E3030, M0009, M3030 and E4002 are GE supplied cables.
  - **All other wiring is customer supplied.**
  - All circuits require Ground wires.

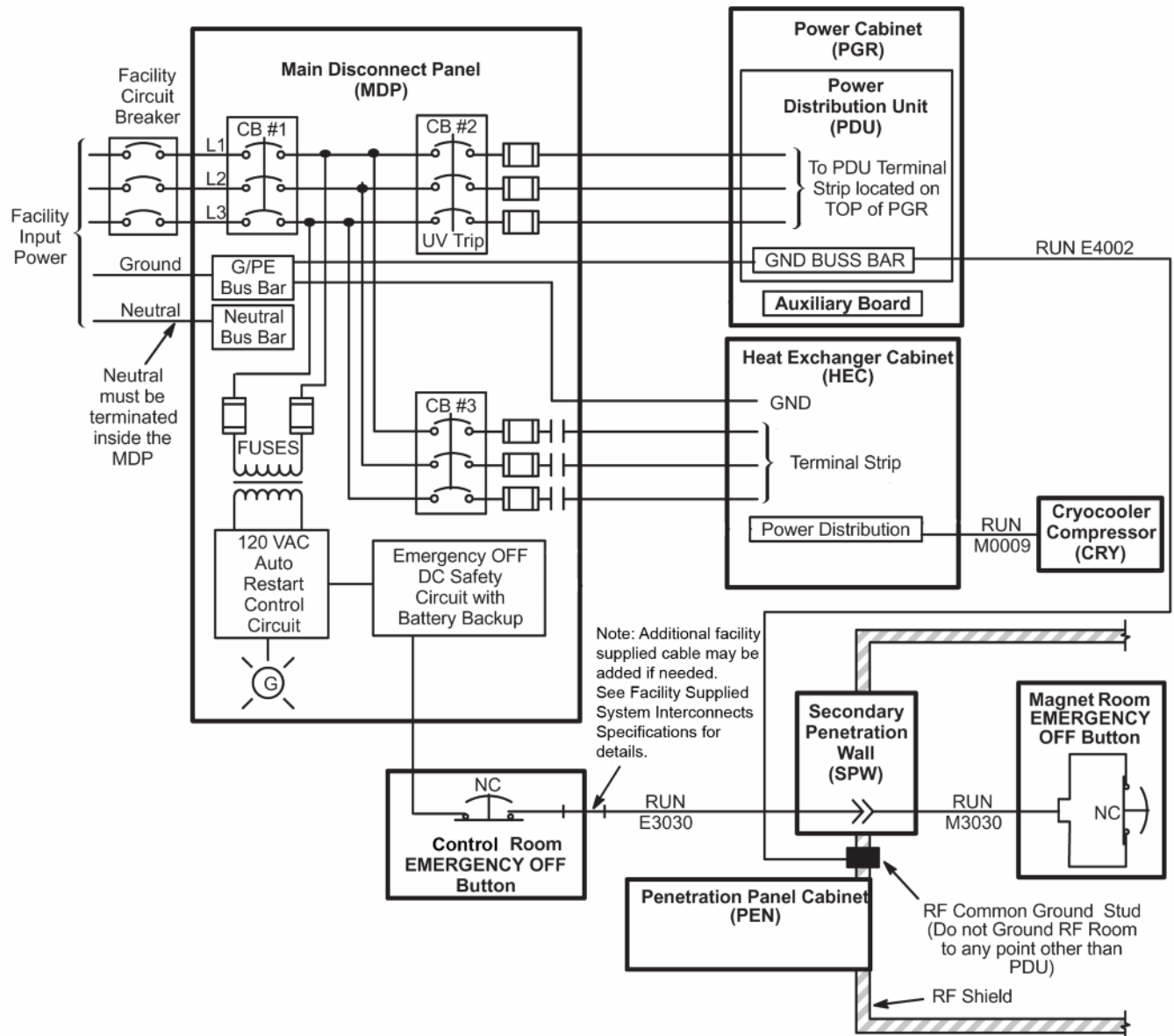

## 12 MR System Shipping and Receiving



### NOTICE

All shipping dimensions and weights are approximate and may vary based on ship-to location, required rigging, or other requirements. Some shipping or access routes may have requirements in addition to those listed in this section. Contact the GE Healthcare Project Manager of Installation (PMI) to verify magnet shipping, rigging, and access.

### 12.1 MR System Shipping Requirements

1. The customer must provide an area for unloading system components from the truck and delivering to the MR suite

**NOTE:** Contact the GE Healthcare Project Manager of Installation (PMI) for a list of experienced rigging companies.

2. At delivery, the customer is responsible for ensuring:
  - a. All floors along the route will support the weight of the magnet (GE Healthcare recommends a structural analysis)
  - b. Doors or other openings are sufficiently wide to allow passage
  - c. Sufficient room is provided for any required dollies or rigging

### 12.2 Facility Delivery Route Requirements

The following table lists the delivery dimensions of system components. The delivery route must be planned to accommodate the dimensions listed.

**Table 2-13: Delivery Route Requirements**

| Component   | Width                                                                       |      | Depth |          | Height |          | Weight   |          | Comments                                                                                                                                                                                                                           |
|-------------|-----------------------------------------------------------------------------|------|-------|----------|--------|----------|----------|----------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|             | in                                                                          | mm   | in    | mm       | in     | mm       | lbs      | kg       |                                                                                                                                                                                                                                    |
| Magnet      | See Magnet shipping dimensions <a href="#">Table 2-14</a>                   |      |       |          |        |          |          |          |                                                                                                                                                                                                                                    |
| Cryogen     | Dimensions vary depending on dewar type used. Verify with cryogen supplier. |      |       |          |        |          |          |          |                                                                                                                                                                                                                                    |
| PGR Cabinet | 78.75                                                                       | 2000 | 34.3  | 872      | 83     | 210<br>8 | 345<br>9 | 111<br>5 | Cabinets are moved with dollies attached to each side (adding 20 in. (520 mm) to the width and 180 lbs (82 kg) to the weight). Cabinets must be raised to remove the pallet but may be lowered almost to floor level while moving. |
| HEC Cabinet | 55.2                                                                        | 1400 | 34.3  | 872      | 75     | 190<br>0 | 113<br>0 | 513      |                                                                                                                                                                                                                                    |
| PEN Cabinet | 44.1                                                                        | 1120 | 39.4  | 100<br>0 | 76     | 193<br>0 | 819      | 371      |                                                                                                                                                                                                                                    |

### 12.3 MR System Component Shipping Specifications

MR system component shipping dimensions and weight are listed below:

Table 2-14: MR System Component Shipping Specifications

| MR Component                                                                       | Approximate W x D x H |                    | Approximate Weight                                                     |                                                                    | Method Of Shipment                            |
|------------------------------------------------------------------------------------|-----------------------|--------------------|------------------------------------------------------------------------|--------------------------------------------------------------------|-----------------------------------------------|
|                                                                                    | in.                   | mm                 | lbs                                                                    | kg                                                                 |                                               |
| Magnet as shipped (with lifting bars)                                              | 93 x 144 x 107        | 2362 x 3658 x 2718 | Domestic:<br><br>11,735<br><br>International:<br>Add 2,200 for Crating | Domestic:<br><br>5322<br><br>International:<br>Add 998 for Crating | Domestic: Tarped. International: Crate/Pallet |
| Magnet Accessory Equipment                                                         | 48 x 48 x 28          | 1219 x 1219 x 711  | 400                                                                    | 182                                                                | Crate                                         |
| Cryocooler Compressor                                                              | 26 x 28 x 42          | 660 x 711 x 1067   | 275                                                                    | 125                                                                | Pallet with box cover                         |
| Rear Pedestal Assembly with Rear Split Bridge Assembly, Low Profile Carriage Cover | 34 x 58 x 48          | 864 x 1473 x 1219  | 310                                                                    | 132                                                                | Box on pallet                                 |
| Enclosure Top                                                                      | 48 x 36 x 36          | 1219 x 914 x 914   | 30                                                                     | 14                                                                 | Box                                           |
| Enclosure Skirts                                                                   | 40 x 24 x 24          | 1016 x 610 x 610   | 30                                                                     | 14                                                                 | Box                                           |
| Patient Table                                                                      | 94 x 29 x 38          | 2377 x 722 x 952   | 474                                                                    | 215                                                                | Pallet                                        |
| Power Gradient RF Cabinet (PGR) (Domestic)                                         | 67.25 x 41 x 88       | 1708 x 1041 x 2235 | 3654                                                                   | 1657                                                               | Pallet                                        |
| Power Gradient RF Cabinet (PGR) (International)                                    | 68.25 x 42 x 92       | 1734 x 1067 x 2337 | 3654                                                                   | 1657                                                               | Pallet                                        |
| Penetration Panel Cabinet (PEN) (Domestic)                                         | 43.75 x 34.25 x 80.75 | 1111 x 870 x 2051  | 710                                                                    | 322                                                                | Pallet                                        |
| Penetration Panel Cabinet (PEN) (International)                                    | 44.75 x 35.25 x 84.75 | 1137 x 895 x 2153  | 770                                                                    | 350                                                                | Pallet                                        |
| Heat Exchanger Cabinet (HEC) (Domestic)                                            | 45 x 44.75 x 80       | 1143 x 1137 x 2032 | 1075                                                                   | 488                                                                | Pallet                                        |
| Heat Exchanger Cabinet (HEC) (International)                                       | 46 x 45.75 x 84       | 1168 x 1162 x 2134 | 1150                                                                   | 522                                                                | Pallet                                        |
| Secondary Penetration Wall                                                         | 21 x 24 x 64          | 533 x 610 x 1626   | 101                                                                    | 45                                                                 | Pallet with cardboard cover                   |
| SPT Phantom Set                                                                    | 34 x 32.5 x 60        | 864 x 826 x 1524   | 350                                                                    | 159                                                                | On cart casters with box cover                |
| Operator Workspace Cabinet                                                         | 24 x 35 x 31          | 600 x 900 x 780    | 243                                                                    | 110                                                                | Wood pallet with cardboard cover              |
| Operator Workspace Display                                                         | 27 x 33 x 27          | 686 x 838 x 686    | 125                                                                    | 57                                                                 | Pallet                                        |
| Operator Workspace equipment                                                       | 32 x 32 x 23          | 813 x 813 x 584    | 100                                                                    | 45                                                                 | Box                                           |
| Operator Workspace Table                                                           | 45 x 54 x 37          | 1143 x 1372 x 940  | 180                                                                    | 82                                                                 | Box                                           |

Table 2-15: MR System Optional Component Shipping Specifications

|                          |              |                 |     |     |                                          |
|--------------------------|--------------|-----------------|-----|-----|------------------------------------------|
| BrainWaveHW Lite Cabinet | 24 x 23 x 72 | 610 x 584 x 320 | 320 | 147 | On cabinet casters, wrapped with plastic |
|--------------------------|--------------|-----------------|-----|-----|------------------------------------------|

|                              |              |                  |     |     |               |
|------------------------------|--------------|------------------|-----|-----|---------------|
| VibroAcoustic Damping Option | 36 x 65 x 12 | 914 x 1651 x 305 | 575 | 261 | Box on pallet |
|------------------------------|--------------|------------------|-----|-----|---------------|


Table 2-16: MR System Component Replacement Shipping Specifications

| Component                                                         | Approximate Weight lbs (kg) | Overall Dimensions W x D x H in. (mm)        | Comments                                                                                                                                             |
|-------------------------------------------------------------------|-----------------------------|----------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------|
| Split Bridge                                                      | 40 (18)                     | 21.5 x 77.3 x 7 (546 x 1969 x 177.8)         |                                                                                                                                                      |
| Replacement RF Body Coil                                          | 155 (70)                    | 30 x 30 x 60 (762 x 762 x 1524)              | Replacement coil is shipped in a protective case. Weight & dimensions are for coil & case.                                                           |
| Replacement XRM Gradient Coil Assembly on a Shipping Cradle/ Cart | 3194 (1449)                 | 39 x 99.84 x 59 (991 x 2536 x 1499) See Note | Initial Gradient Coil Assembly is shipped installed in the Magnet. Shipping/installation cart is used to install replacement coil assembly only.     |
| Gradient Coil Replacement Tool Kit Crate                          | 750 (340)                   | 30 x 86 x 28 (762 x 2184 x 711)              | Gradient Coil Assembly and shipping cart dimensions are with cart in lowest position. Cart can be adjusted to maximum height of 61.88 in. (1572 mm). |

**Note:**

The replacement Gradient Coil Assembly weight is approximately 2205 lbs (1000 kg), the shipping cradle is 132 lbs (60 kg), and the Gradient Coil Assembly shipping/installation cart weighs 855 lbs (389 kg). The coil assembly outside diameter x length dimensions are 35.7 x 57.2 in. (908 x 1452 mm).

Illustration 2-15: Gradient Coil Cart

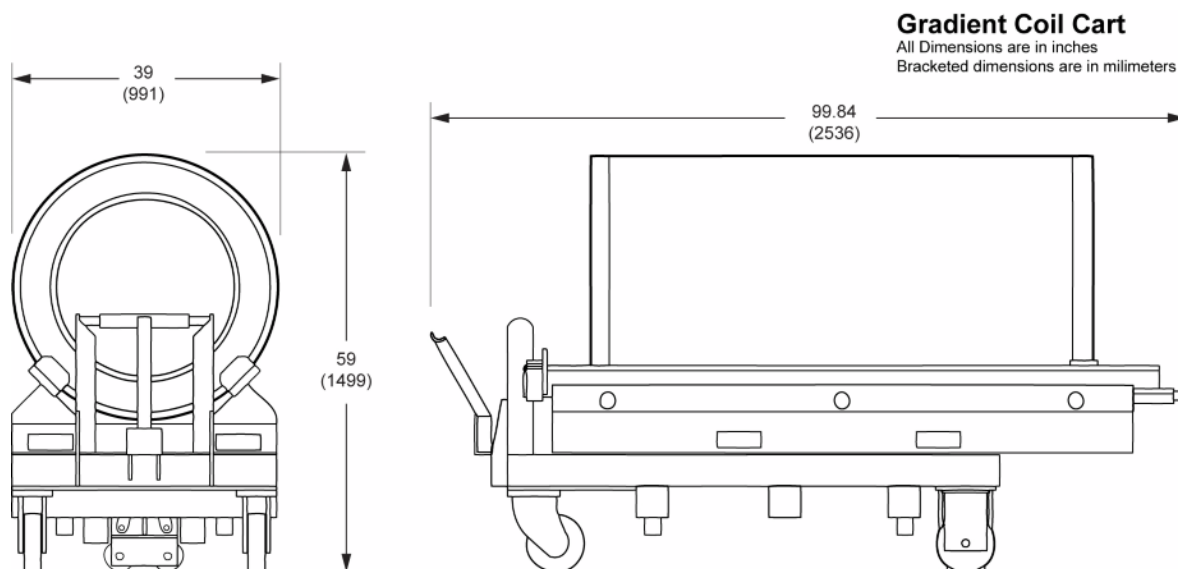




# Chapter 3 Magnet Room

## 1 Magnet Room Introduction

The Magnet Room is best understood as a series of layers, or “rooms within a room.” Each of these rooms has a specific function and associated requirements. All requirements in this chapter must be followed to ensure safe and proper operation of the MR system.

1. The Magnetic shielded room contains the MR Magnet fringe field within a confined space. A site survey is required to determine magnet shield requirements (not all sites require magnetic shielding). Because of the added cost of magnetic shielding, room location should be carefully considered.
2. The Acoustic room is a layer used to help attenuate the noise produced during a scan. An acoustic engineer is strongly recommended to assess the environment.
3. The RF Shielded room is critical to the proper MR system operation. RF shielding prevents interaction of external RF radiation with MR system operation and it also prevents MR system RF radiation from interfering with external systems, such as aircraft control). Special care must be used when installing all fixtures penetrating the RF shield (e.g., vents, electrical conduit, penetration panels) to ensure the integrity of the RF shield is maintained.
4. The Finished room includes the wall coverings, ceiling tile, other fixtures, Magnet (MAG) and Patient Table (PT). When planning the finished room, ensure all building codes (such as maintaining egress routes) are met, items which may cause RF interference, or magnetic items (which could become projectiles when the magnet is installed) are not used.

**Illustration 3-1: Magnet Room Layers**

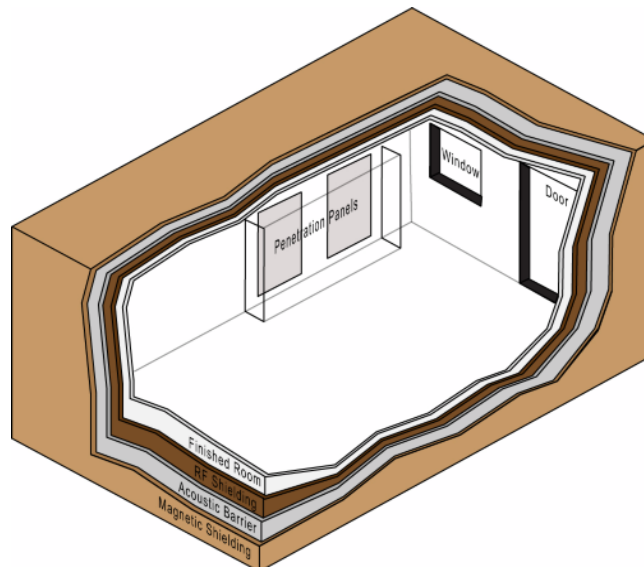

## 2 Magnet Room Structural Requirements

This section lists the structural requirements that must be considered when performing site evaluation and planning of the Magnet room.

### 2.1 Overview

1. When preparing a building plan or evaluating a potential site for an MR system, care should be taken to so the MR suite will not interact with the surrounding environment (i.e., magnetic, acoustic and vibration)
2. The customer is responsible for vibration testing required to verify suitability of a proposed site. All test results and any questions regarding testing, results, or analysis must be forwarded to the GE Healthcare Project Manager of Installation (PMI)

### 2.2 Environmental Steel Limits

A static magnetic field extends in a three-dimensional space around the magnet isocenter. Environmental steel within the static magnetic field affects the uniformity (or homogeneity) of the field. Field uniformity is critical to both image quality and chemical shift analysis (spectroscopy). An analysis of the environmental steel is required in within a 9.84 feet (3 meters) spherical radius of the magnet isocenter. Environmental steel includes pipes, beams, concrete rebar, or any other structural steel in the floors, walls, or ceiling.

The following items must be limited per [Table 3-1](#).

1. Non-movable steel construction material such as wall studs or HVAC components
2. Metallic pipes and drains
3. Steel in the floor in a 10 ft x 10 ft (3.1 m x 3.1 m) area under the Magnet (MAG)

Table 3-1: Steel Mass Limits to Magnet Isocenter


| Limits Of Steel Mass lbs/ft <sup>2</sup> (kg/m <sup>2</sup> ) | Distance From Magnet Isocenter in. (mm) | Distance Below Top Surface Of Floor in. (mm) |
|---------------------------------------------------------------|-----------------------------------------|----------------------------------------------|
| 0 (0)                                                         | 0-45 (0-1143)                           | 0-3 (0-76)                                   |
| 2 (9.8)                                                       | 45-47 (1143-1194)                       | 3-5 (76-127)                                 |
| 3 (14.7)                                                      | 47-52 (1194-1321)                       | 5-10 (127-254)                               |
| 8 (39.2)                                                      | 52-55 (1321-1397)                       | 10-13 (254-330)                              |
| 20 (98.0)                                                     | 55+ (1397+)                             | 13+ (330+)                                   |

## 2.3 Vibration Requirements

Excessive vibration can affect MR image quality. Vibration testing must be performed early in the site planning process to ensure vibration is minimized. Both steady state vibration (exhaust fans, air conditioners, pumps, etc.) and transient vibrations (traffic, pedestrians, door slamming, etc.) must be assessed. Specific requirements for vibration mitigation, include:

1. The Magnet (MAG) cannot be directly isolated from vibration. Any vibration issue must be resolved at the source
2. MR Suite HVAC must have vibration isolation
3. A vibration analysis must be performed at the proposed site with the results (and any mitigation) forwarded to the GE Healthcare Project Manager of Installation (PMI). See the [Chapter 7, MR Site Vibration Test Guidelines](#)
4. A transient vibration test must only be performed after a steady-state test has been performed and all steady-state sources of vibration have been mitigated
5. Transient vibration levels above the specified limits in the [Chapter 7, MR Site Vibration Test Guidelines](#) must be analyzed
6. Any transient vibration that causes vibration to exceed the steady-state level must be mitigated

Illustration 3-2: Magnet Steady State Vibration Specifications

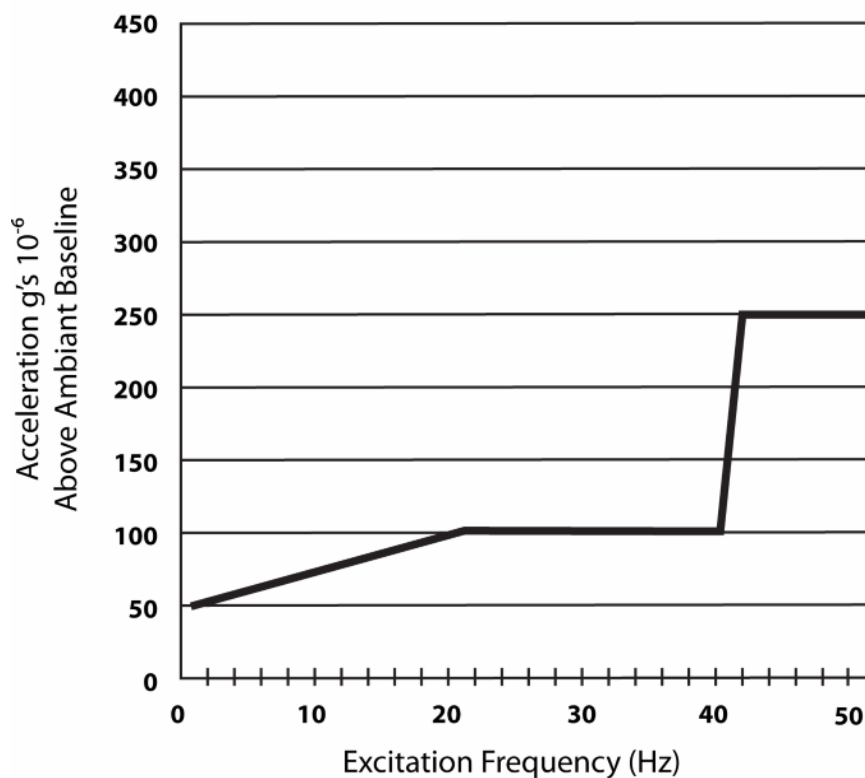Magnet Steady-State  
Vibration Specifications

### 3 Magnetic Shielded Room Requirements



#### NOTICE

All sites, including upgrade sites, must be evaluated for magnetic shielding requirements. Existing magnetic shielding at an upgrade site may not be sufficient for the new system. Contact the GE Healthcare Project Manager of Installation (PMI) to request a site evaluation.

Magnetic shielding prevents interaction between the magnet and nearby sensitive devices. Because of the added cost of magnetic shielding, room location should be carefully considered (not all sites require magnetic shielding). See [Chapter 2, MR Suite Magnetic Field Specifications](#) for detailed magnetic proximity limit information.

1. The GE Healthcare Project Manager of Installation (PMI) works with the customer to coordinate the magnetic shielding site evaluation
2. If required, the GE Healthcare Project Manager of Installation (PMI) coordinates the delivery of the magnetic shielding design
3. The customer is responsible for installation of all magnetic shielding

## 4 Acoustic Room Specifications

The Acoustic room is a layer used to help contain the noise produced during a scan. an acoustic engineer is strongly recommended to assess the environment. This information is only provided as a guide. Actual room noise levels may vary based on room design, optional equipment, and usage.

**Table 3-2: Acoustic Specifications**

|             | Average Sound Pressure Level | Frequency    |
|-------------|------------------------------|--------------|
| Magnet Room | 127 dBA                      | 20 to 20k Hz |

Refer to [Chapter 7, Acoustic Background and Design Guidelines](#) for acoustic design information.

## 5 RF Shielded Room Requirements

The RF Shielded room is critical to the proper MR system operation. RF shielding prevents interaction of external RF radiation with MR system operation (it also prevents MR system RF radiation from interfering with external RF systems, such as aircraft control). Special care must be used when installing all fixtures penetrating the RF shield (e.g., vents, electrical conduit, penetration panels, etc.) to ensure the integrity of the RF shield is maintained.

The RF shielded room can be either a free standing shielded structure or a shielded room within an existing room.

### 5.1 RF Noise Definitions

#### Discrete Interference

Discrete interference is fixed-frequency, narrowband RF noise. Potential sources of discrete interference are radio station transmitters and mobile RF transmitting devices. Magnet room RF shielding prevents external RF energy from entering the room and degrading the MR system RF receivers.

#### Broadband Interference

Broadband interference is caused by electrical discharge within the Magnet room. Potential sources of interference can be reduced by limiting static discharge, ensuring all metal-to-metal contact is tight and secure, and ensuring all electrical and grounding requirements are met.

### 5.2 RF Shield Requirements

#### 5.2.1 Customer Responsibilities

1. The customer is responsible for contracting with an RF Shield vendor to design, install, test, the RF shield (on request, the GE Healthcare Project Manager of Installation (PMI) can supply a list of RF Shield vendors)

**NOTE:** The RF shield may not be in a temperature or humidity controlled environment. Shielding, shield support, and associated components must be installed to prevent degradation over the life of the MR system

2. The customer is responsible for maintenance and repair of RF shielding, to include, but not limited to, RF attenuation, door threshold and door seal, pressure equalization vent operation for the life of the MR System
3. All RF shield designs must be provided to the GE Healthcare Project Manager of Installation (PMI)

#### 5.2.2 Requirements

1. The RF shielding must provide 100dB (up to and including 100MHz  $\pm$  10MHz) planewave attenuation for the entire room including all shield penetration points. Note: 100dB (up to and including 150MHz  $\pm$  10 MHz) is recommended for new construction to accommodate future upgrades
2. The RF Shield vendor installs wave guides including the Cryogen Vent. See [Magnet Room Venting Requirements](#)

3. The customer must work with the RF shield vendor to install any waveguides for optional GE Healthcare (e.g., an MRE air hose) or Non-GE options (functional imaging signal cables)
4. The RF shield must be isolated from ground by more than 1000 ohms
5. The RF shielding must be grounded to the RF common ground stud (which is grounded back to the Power Distribution Unit in the Power, Gradient, RF cabinet)
6. Facility power to the Magnet room must not compromise the RF Shield primary ground
7. RF shield materials must not affect magnetic homogeneity
8. Any moving part (such as doors) must not contain magnetic materials
9. RF shield RF attenuation and ground testing is performed by the customer with a report provided to the GE Healthcare Project Manager of Installation (PMI)

### 5.3 Broadband Interference Avoidance Requirements

1. Ensure all joints and mechanical connections remain secure:
  - a. All solder joints must be clean and properly prepared
  - b. All mechanical fasteners must be sufficiently tightened and secured
  - c. Do not use rivets or self-tapping screws (as these tend to loosen over time due to vibration)
2. Prevent RF shield corrosion:
  - a. Avoid contact between dissimilar metals
  - b. Ensure all joints and seams are properly dressed using proper materials
  - c. Use sacrificial anodes
3. Doors and door frames must be structurally stiff to prevent physical changes to the RF shield

### 5.4 Ambient Radio Frequency Interference (RFI) Requirements

The MR System uses spatially encoded radio frequency information to create the MR image. All sites require a 100 dB RF Shield for the entire frequency range defined in the table below. Listed in the table below are the recommended centerband and bandwidth frequencies to be used when measuring radio frequency interference. This table includes frequency bands important for both imaging and spectroscopy.

**Table 3-3: Radio Frequency Survey Specifications**

| Isotope          | Bandcenter MHz | Bandwidth Hz |
|------------------|----------------|--------------|
| <sup>1</sup> H   | 63.86          | 916,138      |
| <sup>19</sup> F  | 60.12          | 981,882      |
| <sup>31</sup> P  | 25.88          | 390,296      |
| <sup>23</sup> Na | 16.90          | 242,773      |
| <sup>13</sup> C  | 16.06          | 233,925      |

RFI site surveys are to be performed by cycling through the preceding frequency bands and a broad band range up to 100MHz  $\pm$  10MHz (up to 150MHz  $\pm$  10MHz is recommended for new sites to accommodate upgrades). Special emphasis, however, should be placed on the 1H band since



this is used in proton imaging. The RFI site survey should be performed for a length of time necessary to determine, within a reasonable degree of certainty, that the RFI noise at the site will not exceed the 100 dB attenuation provided by the RF shielded room. Note that any RFI site survey no matter how thorough, will not preclude the possibility of future or unmeasured RFI caused by new or intermittent sources.

The ambient RF noise measured should be less than 100 millivolt per meter (100 dB microvolt per meter). When a RFI site survey is required, it must be completed before the purchase and installation of the RF shielded room.

To ensure that 100 millivolt (or greater) RF noise peaks outside the bandwidths specified above do not actually extend into these bandwidths and exceed the 100 millivolt limit, adjust the resolution of the test equipment (spectrum analyzer) according to the equation:

$$BW \text{ (resolution)} = f_0 / 50$$

where: BW = Bandwidth (resolution)

$f_0$  = Center frequency (for  $^1\text{H}$ : at 1.5 Tesla 63.86 MHz)

## 5.5 Magnet Mounting Requirements


The following requirements are applicable only if the VibroAcoustic dampening option is not included with the MR system.

1. The RF Shield vendor must design and install magnet anchor bolts
2. All four Magnet feet must be anchored to the floor (refer to [Illustration 3-3](#))
3. Magnet anchors must be installed before the Magnet is delivered
4. Refer to [RF Shielded Room Requirements](#) for Magnet Feet mounting hole locations
5. Magnet anchors must not contact floor rebar or other structural steel
6. Magnet anchors must electrically contact the RF shield at point of entry
7. Magnet anchors must have the following properties:
  - a. Anchors must be two-part assembly (male/female)
  - b. Female side must be expansion- or epoxy-type
  - c. Male side must be a bolt or threaded rod with appropriate-sized nut (bolt or rod must be removable--not epoxied or cemented in place)
  - d. Anchors must be electrically conductive
  - e. Anchors must be non-magnetic
  - f. Anchors must not induce galvanic corrosion with the RF shield
  - g. Anchors must be commercially procured
  - h. Anchor rods diameter must be between 0.625 in (M16) and 1.25 in (M32)
  - i. Anchors must meet the following clamping force: 2,500  $\pm$  200 lbs (11,100  $\pm$  900 N)

**NOTE:** Refer to [Chapter 7, Selecting Magnet Anchor Size](#) for an example of how to properly size the Magnet anchor

8. The RF shield vendor must perform the following tests before installation. Results must be provided to the GE Healthcare Project Manager of Installation (PMI)
  - a. A pull test on each anchor (equal to the clamping force)
  - b. RF integrity test
  - c. Ground impedance test

**Illustration 3-3: Magnet or Dock Anchor Mounting Details**

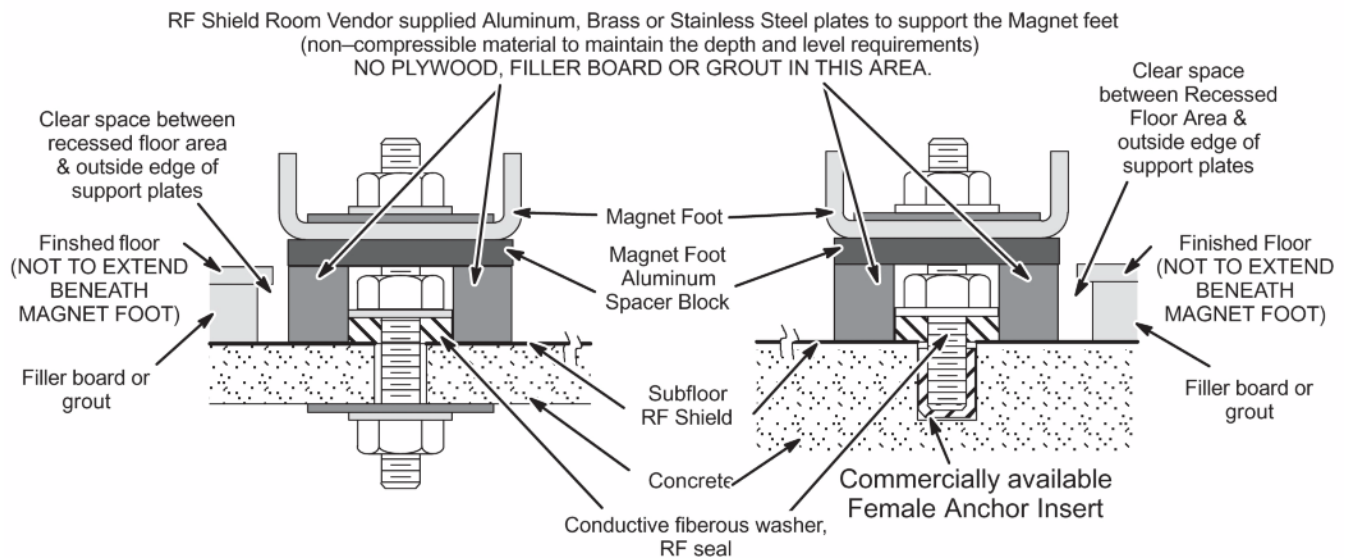

**Note:** For sites with RF Shield on top of subfloor, the RF Shield needs to be recessed to the concrete level to provide a proper RF Seal.

## 5.6 Dock Anchor Mounting Requirements

1. The RF Shield vendor must design and install magnet anchor bolts
2. The dock anchor hole must be drilled after the Magnet is installed
3. The dock anchor must not contact floor rebar or other structural steel
4. The dock anchor must electrically contact the RF shield at point of entry
5. The dock anchors must have the following properties (refer to [Illustration 3-3](#)):
  - a. Anchors must be two-part assembly (male/female)
  - b. Female side must be expansion- or epoxy-type
  - c. Male side must be a bolt or threaded rod with appropriate-sized nut (bolt or rod must be removable--not epoxied or cemented in place)
  - d. Anchors must be electrically conductive
  - e. Anchors must be non-magnetic
  - f. Anchors must not induce galvanic corrosion with the RF shield

- g. Anchors must be commercially procured
  - h. The anchor rod hole clearance in the dock anchor base is 0.43 in (11). The anchor rod diameter must be sized appropriately
  - i. Anchors must meet the following clamping force:  $600 \pm 100$  lbs ( $2,700 \pm 450$  N)
6. The RF shield vendor must perform the following tests. Results must be provided to the GE Healthcare Project Manager of Installation (PMI)
- a. A pull test on each anchor (equal to the clamping force)
  - b. RF integrity test
  - c. Ground impedance test

## 6 Finished Room Requirements



**Personnel Injury or Equipment Damage**

**Metallic objects may become projectiles if not properly secured.**

**Remove or properly secure any metallic objects within the finished room.**

### 6.1 Walls

1. GE Healthcare recommends wall coverings to protect the RF shielding (wall coverings may be removable, if required, to access the RF shield)
2. Metallic electrical conduit inside walls and ceilings may be used. Conduit for receptacles must be metallic
3. An enclosure (i.e., PEN closet) must be provided to restrict access to the PEN wall and for storage of excess interconnections
  - a. The PEN closet must have a mechanical locking mechanism to restrict access to the PEN panels
  - b. The PEN closet must enclose the minimum service area in the Magnet room as shown in [Chapter 4, PEN and SPW Wall Opening Requirements](#).
  - c. The PEN closet may be expanded to provide an area for excess cable storage with the following requirements:
    - i. Excess cable must not be stored within the minimum service area
    - ii. Excess cable must not interfere with access or servicing of the PEN panel or SPW
    - iii. The area within the PEN closet to store the cable should be sized to accept a 22 in (559 mm) cable loop (2x the minimum bend radius of the largest cable)

### 6.2 Doors and Magnet Access Openings


1. The finished opening of the main door must be at least 43 in. (1092 mm) to allow for helium dewars and patient tables
2. The finished opening of the main door must be positioned to allow the Gradient Cart to move into the room
3. Threshold height must not exceed 1 in. (25 mm) on both sides of the door with a maximum 10 degree threshold inclination
4. Patient viewing windows recommended dimensions are 48 in. wide by 42 in. high (1219 mm x 1067 mm) and 72 in. (1829 mm) above the finished floor

**NOTE:** IEC requires the patient, while in the bore, be in full view of the operator. GE Healthcare recommends using a window, although other means (e.g., camera and video display) may be used as long as all IEC requirements are met
5. The magnet delivery requires opening into the room to allow access for the magnet and any rigging

### 6.3 Finished Ceiling

1. The customer is responsible for the finished ceiling, including any cable enclosures
2. The illustration below shows the maximum dimensions of a customer supplied enclosure at the magnet opening (the ceiling opening size can vary based on customer requirements)

### Illustration 3-4: Ceiling Interface

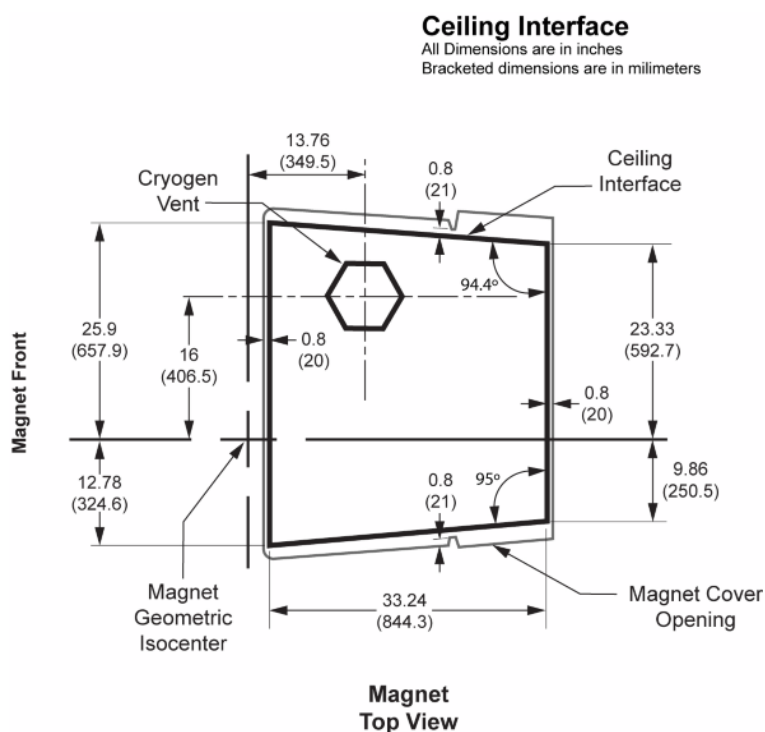


**NOTE:** The illustration above is provided for reference only. Actual dimensions may vary due to construction tolerance limitations and magnet placement tolerance limitations.

## 6.4 Magnet Room Floors


1. The finished floor must support the weight of all components (e.g., patient table, gradient coil replacement cart) throughout operation and service life
2. The finished floor must be water resistant to protect the subfloor and shielding from water damage
3. The customer is responsible for providing flooring to meet the maximum 8 kV ESD resistance requirement for GE Healthcare equipment
4. Magnet, Enclosure, and Patient Table areas (shown below) must be level to 0.3125 in. (8 mm) between high and low spots over any 120 in. (3048 mm) distance
5. If the Magnet is mounted directly to the floor, see [RF Shielded Room Requirements](#) for mounting bolt requirements
6. When using the GE Healthcare supplied VibroAcoustic Dampening option, use 24 in. x 65 in. (610 mm x 1651 mm) for each pad to calculate floor loading

7. RF shield seams, joints, or overlaps must not be located under the VibroAcoustic mats

**Illustration 3-5: Magnet Room Floor Levelness Area**

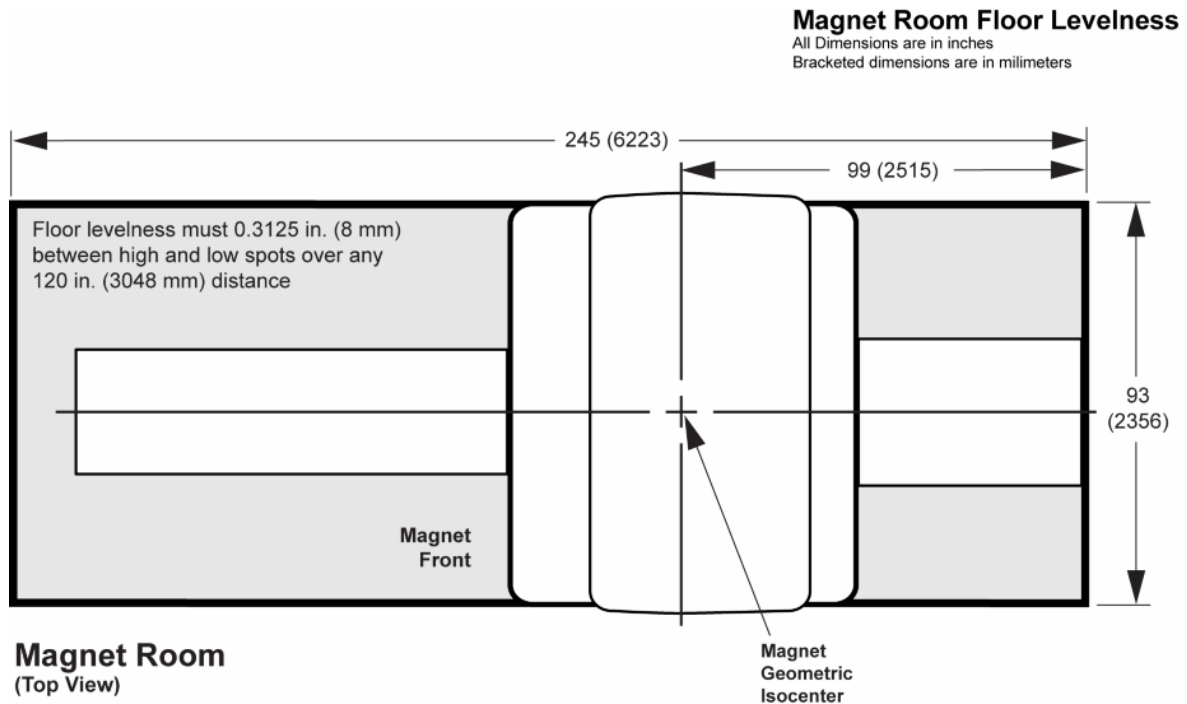

**Illustration 3-6: Magnet Mounting Detail**

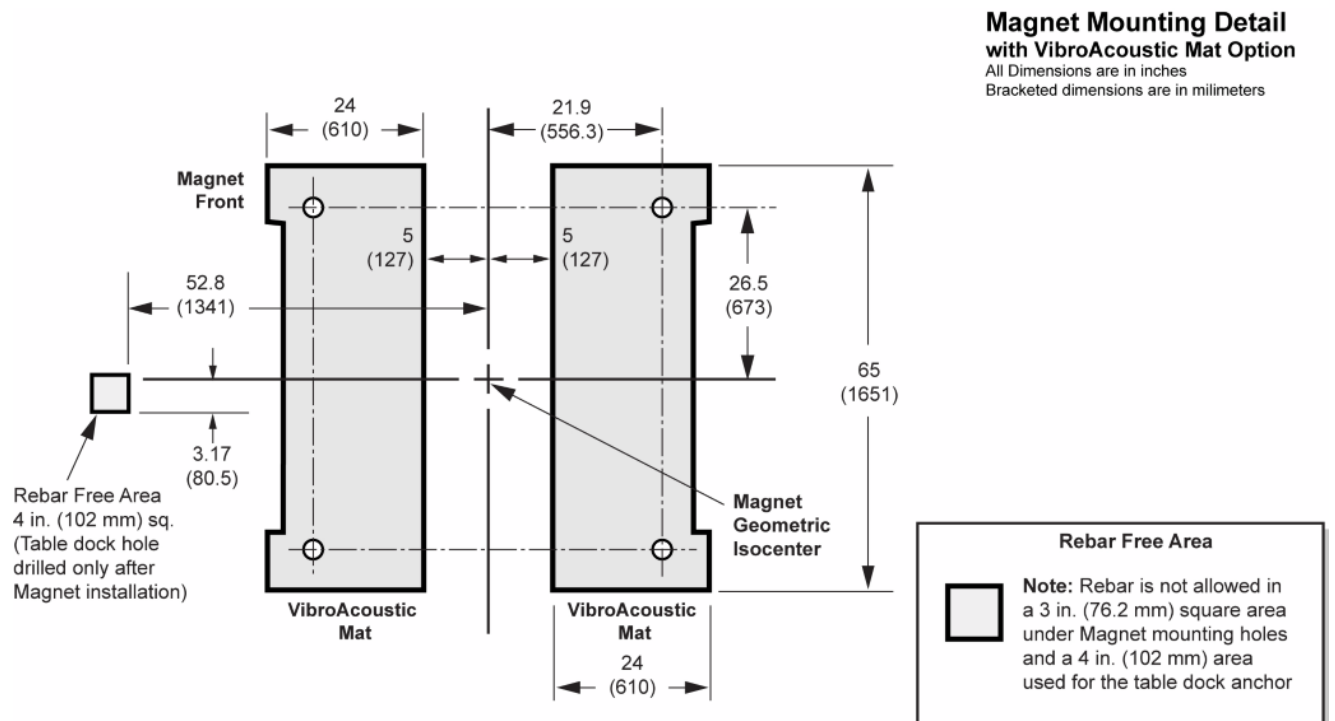


## 7 Magnet Room Equipment Specifications

### 7.1 Magnet (MAG) Specifications

Magnet weight (installed): 11,475 lbs (5,205 kg)

Rear Pedestal weight: 212 lbs (96 kg)

Illustration 3-7: Magnet (MAG) and Rear Pedestal

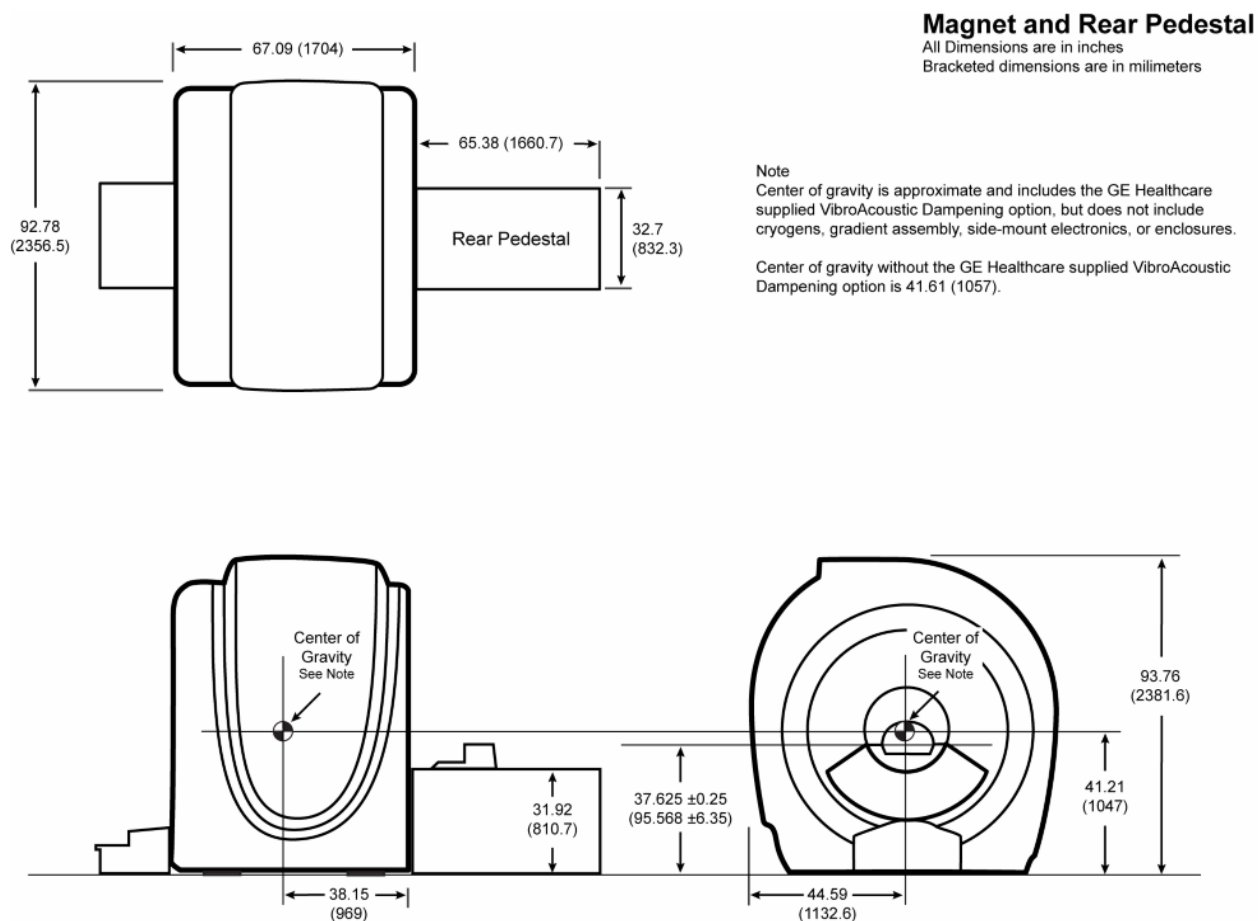


### 7.2 Patient Table (PT) Specifications

Weight, empty: 350 lbs (159 kg)

Weight with maximum patient weight of 500 lbs (227 kg): 850 lbs (386 kg)

Illustration 3-8: Patient Table (PT)

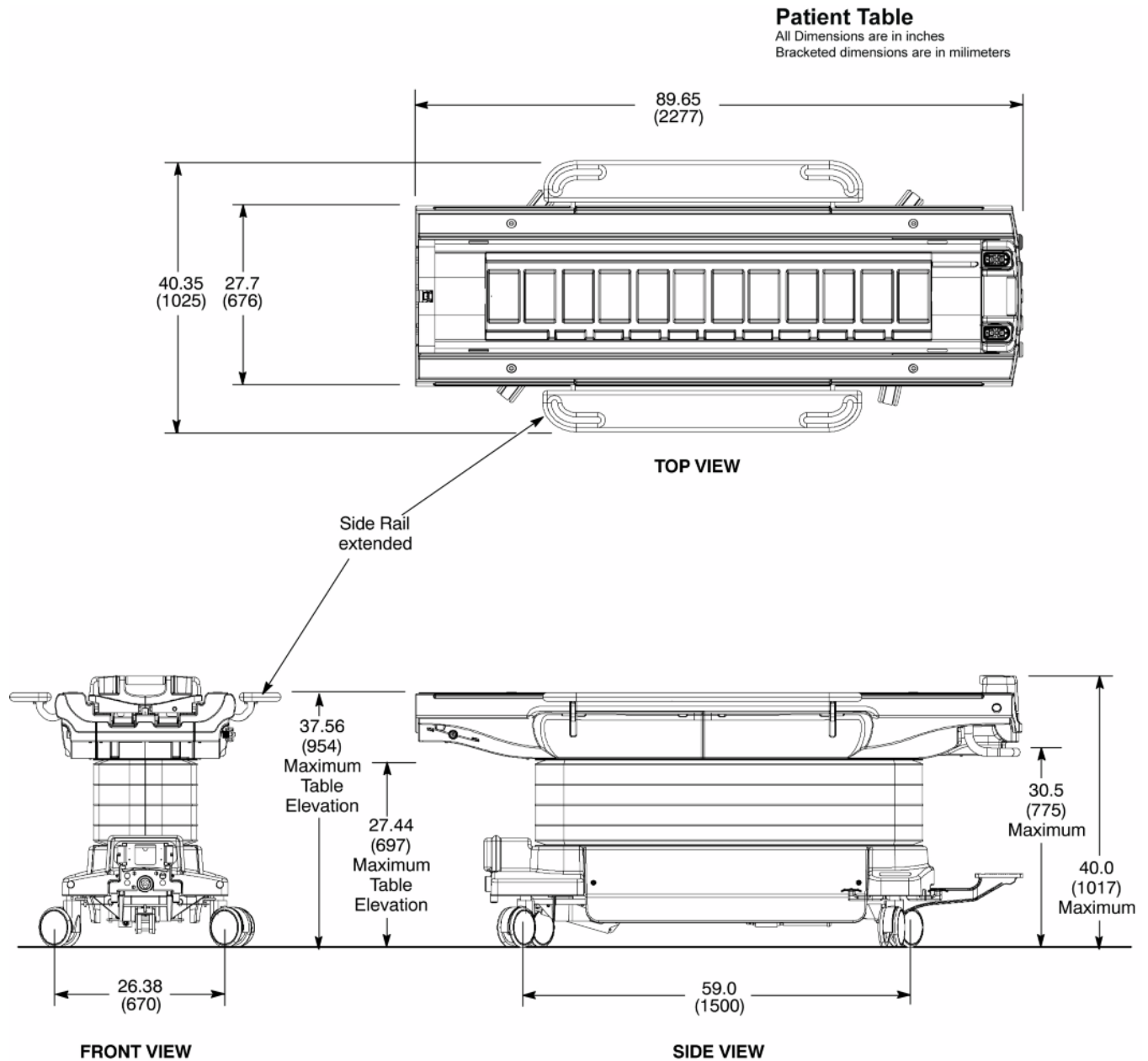

### 7.3 Magnet Rundown Unit (MRU) Specifications and Requirements

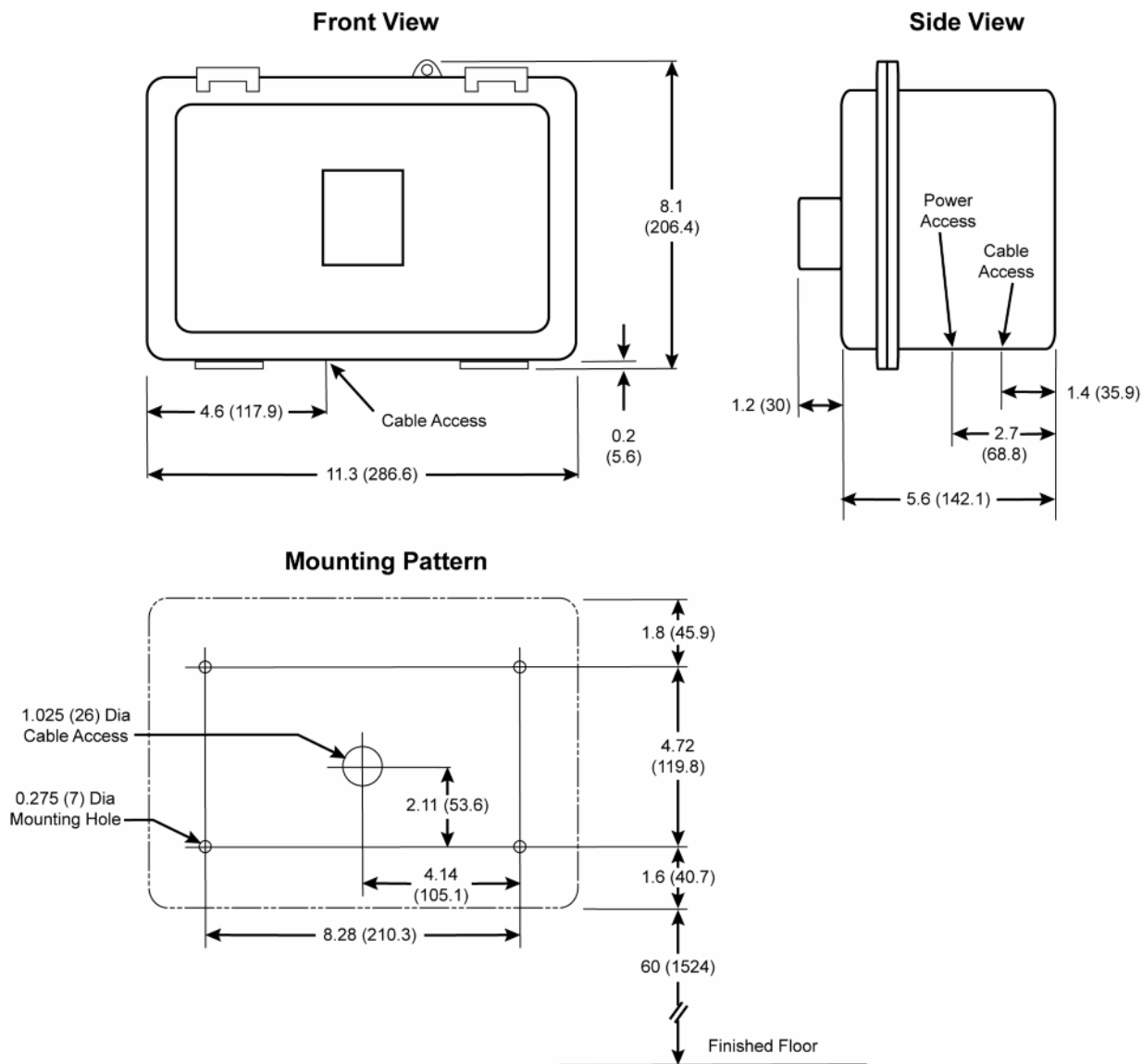
1. Location: The bottom edge of the MRU must be mounted 60 in. (1524 mm) above the Magnet room floor near the front of the magnet enclosure
2. Weight: 7 lbs (3.2 kg)
3. Magnetic field limit: 200 gauss (20 mT)
4. The MRU is installed by the facility contractor
5. The MRU requires the following facility supplied power:



Table 3-4: MRU Facility Power Requirements

| Parameter           | Requirements                                                                                  |          |
|---------------------|-----------------------------------------------------------------------------------------------|----------|
| Voltage / Frequency | 100-120 VAC                                                                                   | 50/60 Hz |
|                     | 200-220 VAC                                                                                   | 50/60 Hz |
| Phase               | 1                                                                                             |          |
| Maximum Amps        | 1.0                                                                                           |          |
| Connection type     | Hard wired in unit. Plastic conduit is required (1/2 in. PVC Schedule 40 Conduit recommended) |          |
| Availability        | Continuous                                                                                    |          |
| Circuit Breaker     | Dedicated AC disconnect required for both live and neutral connections                        |          |

Illustration 3-9: Magnet Rundown Unit (MRU)

**Magnet Rundown Unit (MRU)**All Dimensions are in inches  
Bracketed dimensions are in millimeters

## 7.4 Oxygen Monitor Sensor Specifications

See [Chapter 5, Oxygen Monitor \(OXY\)](#)

## 8 Magnet Room Venting Requirements

### 8.1 Overall Requirements

1. The Magnet room must have four independent venting systems:
  - a. HVAC vents to maintain airflow and proper temperature in the Magnet room
  - b. An exhaust vent system to remove any cryogenic gasses released into the room
  - c. A pressure equalization vent
  - d. A cryogenic vent to allow helium gas to escape during a quench
2. All vent components must be accessible for inspection, cleaning, or maintenance

### 8.2 HVAC Vent Requirements

1. All components of the HVAC system are supplied by the customer
2. RF Shield vendor must install open pipe or honeycomb HVAC waveguides
3. Waveguides must be nonmagnetic and electrically isolated
4. Incoming air must contain at least 5% air from outside the Magnet room (refer to [Illustration 3-10](#) and [Illustration 3-11](#))

### 8.3 Emergency Exhaust Vent Requirements

1. All components of the exhaust vent system are supplied by the customer
2. The exhaust vent system must be operational before the magnet is installed
3. The exhaust intake vent must be located near the magnet cryogenic vent at the highest point on the ceiling
4. The Magnet room exhaust fan and exhaust intake vent must have a capacity of at least 1200 CFM (34 m<sup>3</sup>/minute) with a minimum 12 room air exchanges per hour
5. The exhaust fan must be placed above RF shielding located outside 10 gauss (1mT) and with appropriate waveguide
6. The system must have a manual exhaust fan switch near the Operator Workspace (OW) and in the Magnet room near the door (the switches must be connected in parallel)
7. All system components must be accessible for customer inspection, cleaning, and maintenance

Illustration 3-10: Magnet Room Exhaust Fan Schematic

**NOTE:**

All items shown are supplied and installed by Customer or Contractors.


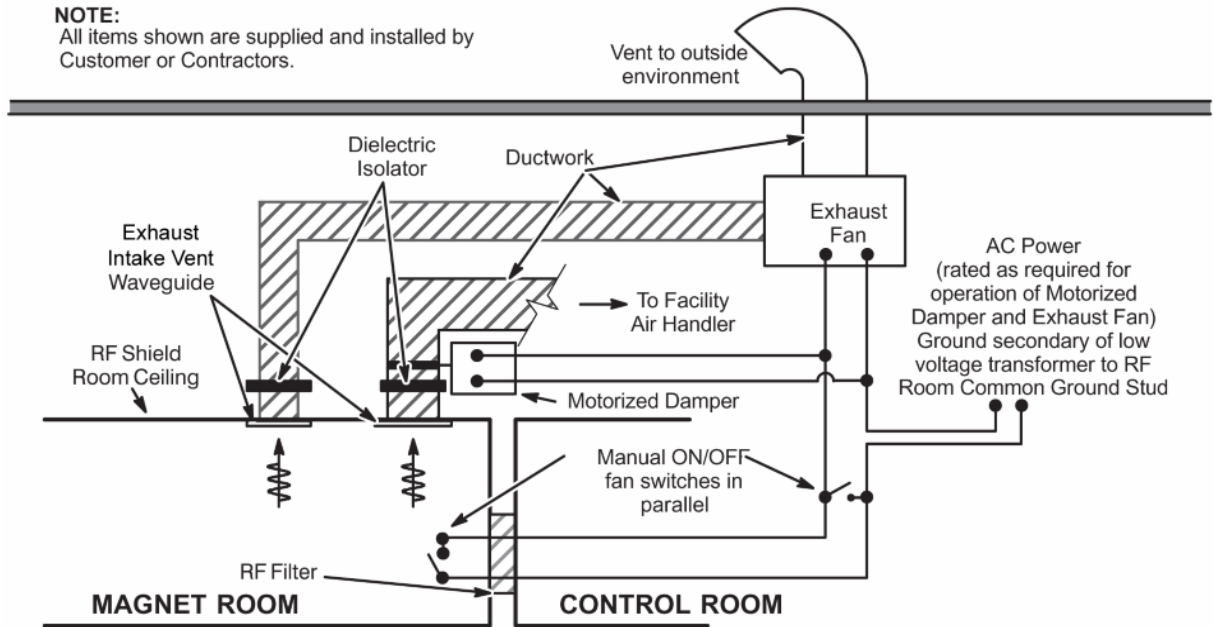


Illustration 3-11: Magnet Room Exhaust Fan Schematic with Optional Oxygen Monitor (OM)

**NOTE:**

GE Option supplied and installed equipment. All other items supplied and installed by Customer or Contractors.

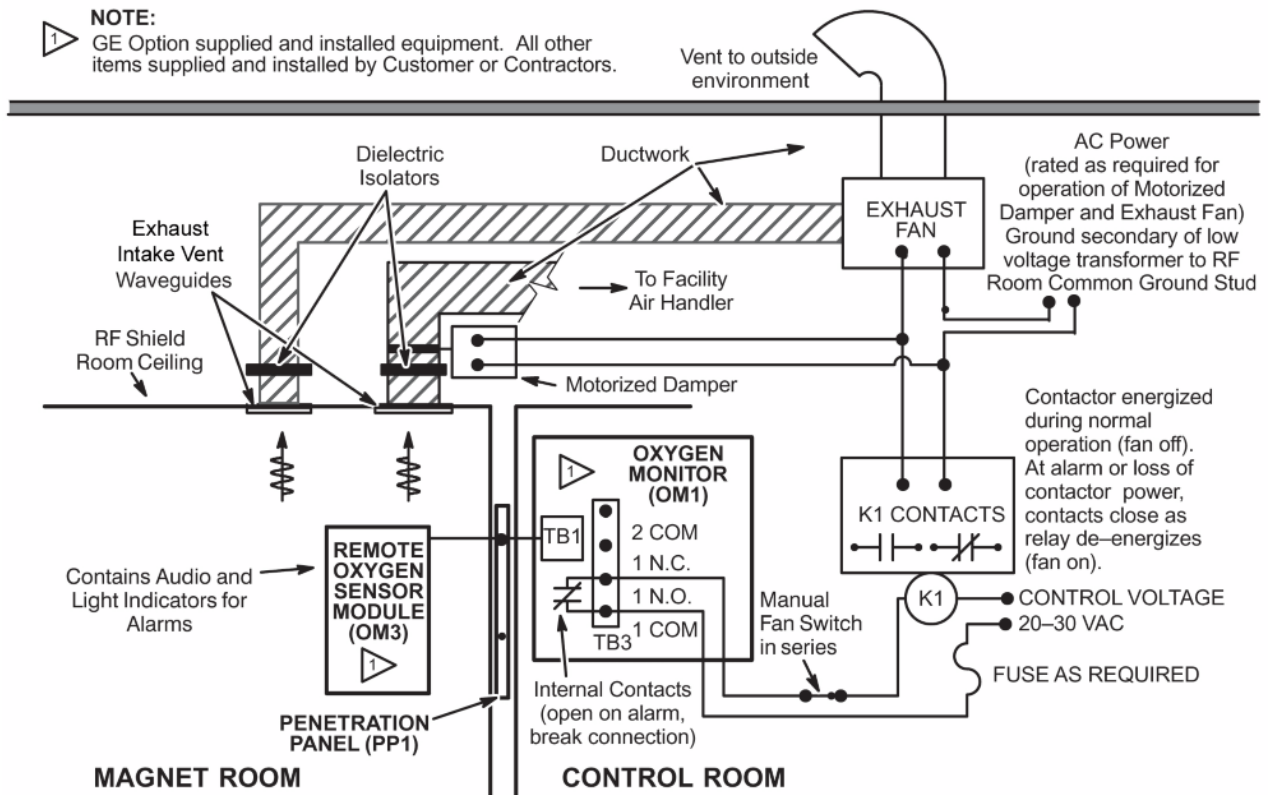

## 8.4 Pressure Equalization Vent Requirement

An additional pressure equalizing waveguide vent in the magnet room ceiling or in the wall (with waveguide top edge located at ceiling) to prevent air pressure from restricting movement of the Magnet room door. The vent minimum size must be either:

1. 24 in. x 24 in. (610 mm x 610 mm)
2. 576 sq in. (372,100 sq mm)

## 8.5 Cryogenic Venting

The MR System requires an outside vent to allow helium gas to escape during magnet quench. The vent must be properly installed to ensure the gas does not escape into the building. The vent must be connected to the magnet within 24 hours of installation in the magnet room.

Note the following:

1. All pipe or tube dimensions specified in this document are outside diameters unless otherwise noted
2. Some requirements in this document may not be applicable to Mobile, Transportable, or Relocatable configurations
3. See [Chapter 7, Magnet Cryogenic Venting Pressure Drop Reference Tables](#) to calculate pressure drop for a specific system

Table 3-5: Magnet Cryogen Specifications

| Magnet Types    | Helium Volume gallons (liters) | Peak Helium Flow During Quench ft <sup>3</sup> per min (m <sup>3</sup> per min) | Magnet Vent Pipe OD inches (mm) |
|-----------------|--------------------------------|---------------------------------------------------------------------------------|---------------------------------|
| DVw 1.5T Magnet | 468 (1770)                     | 2737 (77.5) [Gas]                                                               | 8 (203.2)                       |

### 8.5.1 Interior Vent Requirements

#### 8.5.1.1 General

1. Do not remove or modify the vent adaptor bolted to the magnet
2. The GE Healthcare supplied vent tube must be bolted directly to the magnet vent adaptor bolt flange
3. The vent tube may be cut to a minimum of 4 inches (10.6 cm) from the top of the flange
4. The vent must be located within 0.25 in (6.35 mm) of the location (in relation to isocenter) shown in [Illustration 3-13](#)

#### 8.5.1.2 Vent Size

The total pressure drop of the cryogenic vent system (from the magnet vent interface to the termination point outside the building) must be less than 17 psi (117.2 kPa). The pressure drop of the RF shield waveguide must be included in the overall calculation.

Refer to Cryogen Vent System Typical Information for additional calculation specifications.

### 8.5.1.3 Vent Materials

1. The vent material must be one of the following materials with the wall thickness indicated:
  - a. SS 304: Minimum 0.035 in. (0.89 mm); Maximum 0.125 in. (3.18 mm)
  - b. AL 6061-T6: Minimum 0.083 in. (2.11 mm); Maximum 0.125 in. (3.18 mm)
  - c. CU DWV, M or L: Minimum 0.083 in. (2.11 mm); Maximum 0.140 in. (3.56 mm)
2. Either tubes or pipes may be used and must be seamless or have welded seams
3. Corrugated pipe or spiral duct must not be used
4. A bellows pipe less than 1 ft (30 cm) in length may be used as a thermal expansion joint
5. The vent pipe must withstand a maximum pressure of 35 psi (241.4 kPa)
6. Waveguide vent material must match the outside diameter of the magnet vent

### 8.5.1.4 Support

1. The vent support assemblies must support the entire vent system and 1850 lbs (8229 N) helium flow reaction force at vent elbows
2. Any vent support between the magnet interface and the isolation joint at the waveguide must be electrically isolated to ensure the integrity of the RF shield
3. The Ventglas joint must not be used as a vent system support

### 8.5.1.5 Construction

1. A single dielectric break in the vent system (i.e., Ventglas) is required in the Magnet room to ensure the integrity of the RF shield for the Magnet Room
  - a. Ensure the gap between the pipes is 1.0  $\pm$  0.25 inch (25.4  $\pm$  6 mm)
  - b. The outside diameter of the waveguide must match the outside diameter of the GE vent tube within  $\pm$  0.125 in. (3 mm)
  - c. The Ventglas joint must be accessible for annual inspection or maintenance
  - d. Do not use rubber soil pipe
2. The Ventglas joint may also serve as a thermal expansion joint
3. All joints except the single Ventglas joint must be welded or brazed
4. All isolation/thermal expansion joints (except the Ventglas joint) must be rated to 4.5 K (-451°F or -268°C) and 35 psi (241.4 kPa)
5. The vent system must be insulated with 1.5 inch (38 mm) thick flexible unicellular insulation to prevent condensation during magnet ramping. Exposed insulation must be covered with a white PVC jacket
6. The installation of the vent pipe section between the waveguide and the supplied magnet vent pipe must not result in any vertical or side load to the waveguide and the magnet vent pipe (e.g., Lorenz clamps must not be used because the rigid connection will transfer load to the magnet vent pipe)

## 8.5.2 Exterior Vent Requirements

### 8.5.2.1 Vent Support

1. The vent support assemblies must be able to support the entire vent system and 1850 lbs (8229 N) helium flow reaction force at vent elbows
2. The Ventglas joint must not be used as a vent system support

### 8.5.2.2 Vent Construction

1. The vent must be routed as directly as possible to the outside. If used, elbows must be standard or long sweep
2. Expansion/contraction joints must be provided for temperature decrease from ambient to 4.5 K (-451°F or -268 °C)
3. A Ventglas joint must be installed at the waveguide. The joint gap must be  $1.0 \pm 0.25$  inch ( $25 \pm 6$  mm)
4. The GE Healthcare supplied Ventglas and clamps are for 8 in. (203 mm) pipe only
5. The Ventglas joint must be accessible for inspection or maintenance
6. All components must be rated to withstand the helium flow reaction force at temperatures from ambient to 4.5 K (-451°F or -268°C)
7. Electromechanical fire dampers must not be used. Fusible link fire dampers may be used (with annual inspection)
8. Exit of vent must be covered to prevent ingress of weather elements (i.e., rain, snow, hail, sand, etc.) and foreign material debris (i.e., leaves, bird nests, etc.)
9. Condensate must be prevented from pooling inside any section of the venting system

### 8.5.2.3 Vent Exit



#### **CAUTION**

##### **Cryogenic Burns or Asphyxiation**

During a quench, extremely cold gas or particles are released from the cryogenic venting system. A Quench may occur at any time.

**Ensure access to cryogen vent exhaust area is restricted and the released gas does not reenter the building. Refer to the specifications below.**

1. An exhaust area in front of the vent 20 feet (6.1 m) long by 15 feet (4.6 meters) wide:
  - a. Must be restricted and labeled with appropriate warning signs
  - b. Must not include air intake vents
  - c. Must not include any personnel, building components, or objects (movable or stationary)
2. For a rooftop exit:
  - a. Use either a horizontal exhaust vent with a 90° elbow and minimal pressure drop or other low pressure drop, high flow rate roof cap

- b. The bottom of the 90° elbow must be at least 3 feet (0.9 meters) above the roof deck (or higher if at risk of being blocked by drifting snow, sand, etc.)
  - c. The outlet must be covered with a 0.5 inch (12.7 mm) square screen mesh
  - d. The exhaust vent must be included in the pressure drop calculation
3. For a sidewall exit:
- a. Use an exhaust vent with a 45° elbow (with a deflector rated for the helium reaction force) and no restriction in gas flow
  - b. The exhaust exit must be at least 12 feet (3.66 meters) above the ground
  - c. The outlet must be covered with a 0.5 inch (12.7 mm) square screen mesh
  - d. The outlet must be covered to prevent foreign material from entering or blocking the opening (e.g., louvers, etc.)
  - e. The exhaust vent must be included in the pressure drop calculation



## NOTICE

The minimum ceiling height dimensions shown below are for vent installation only. Overall ceiling height must also take into account minimum service area requirements [Chapter 6, MR Suite Minimum Room Size Requirements](#) and Interconnect routing requirements [Chapter 2, MR System Interconnects Routing Requirements](#).



Illustration 3-12: Cryogenic Vent Installation Detail


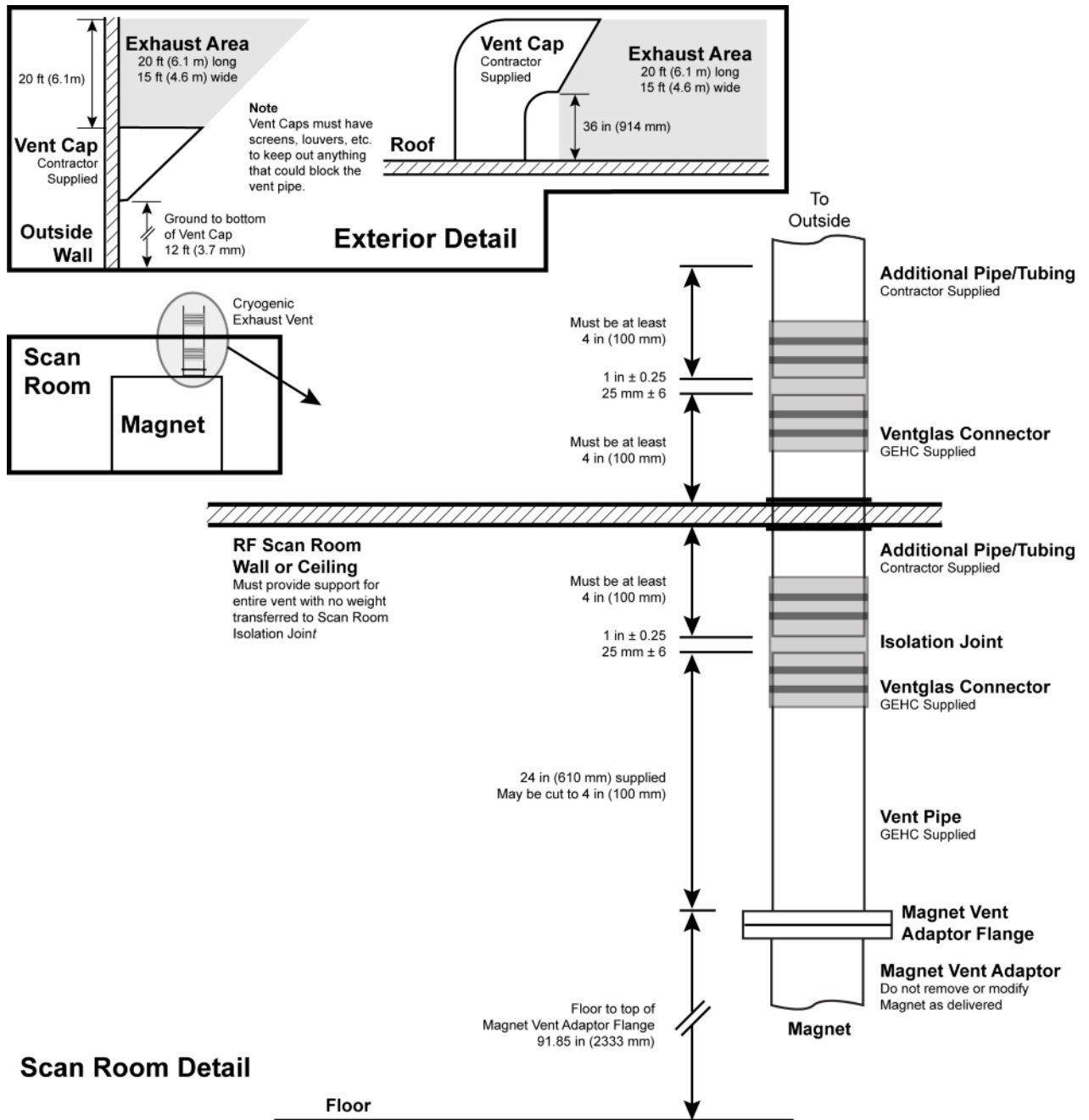


Illustration 3-13: Magnet Cryogenic Vent Location

**Magnet Cryogenic Vent Location**

All Dimensions are in inches  
Bracketed dimensions are in millimeters

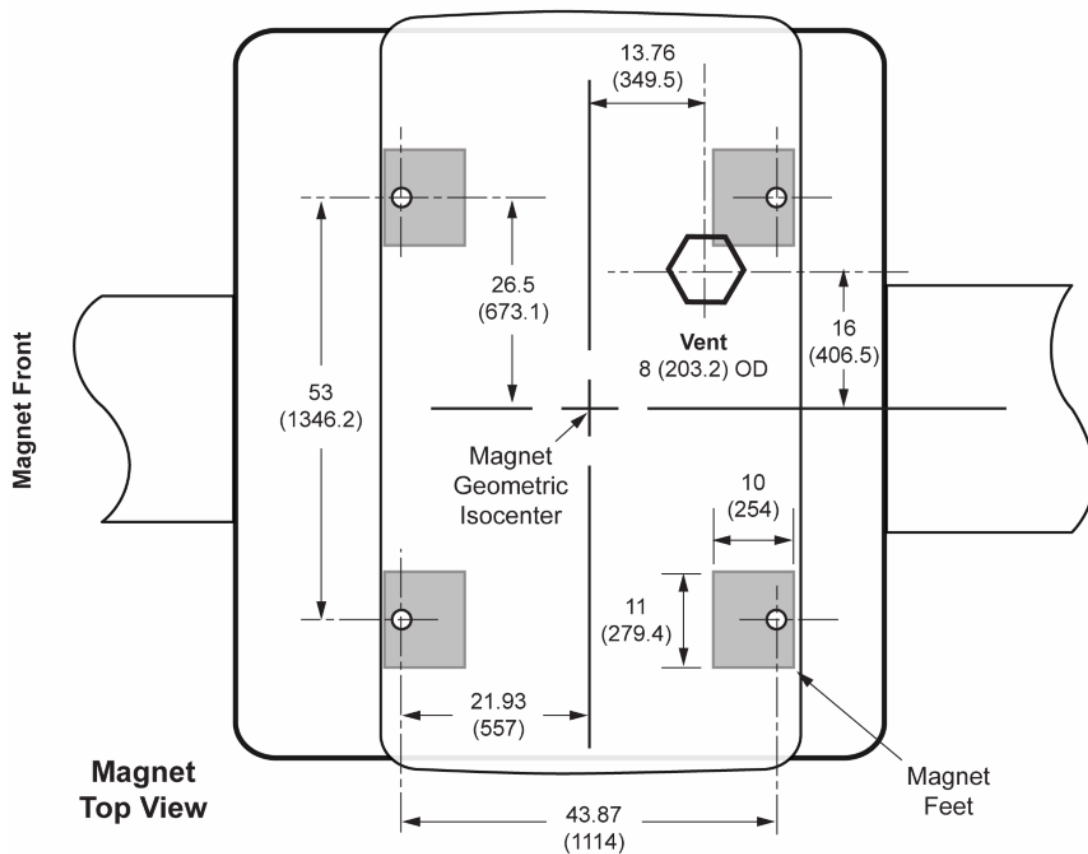

## 9 Magnet Room Electrical and Grounding Requirements

### 9.1 Electrical Line and Filter Requirements

1. RF Shielded Room vendor and electrical contractor must coordinate installation of all electrical lines through the RF shielding
2. The RF Shielded Room vendor must supply electrical line filters for all lines through the RF shielding (excluding electrical lines through the GE supplied Penetration panels) to ensure compliance with the RF Shielded Room attenuation requirements
3. Electrical line filters must be located outside the 200 gauss (20 mT) line


### 9.2 Lighting Requirements

1. All lighting fixtures and associated components must meet all RF Shielded Room and RF Grounding requirements (e.g., track lighting is not recommended due to possible RF noise)
2. All lighting must use direct current (the DC must have less than 5% ripple)
3. 300 lux must be provided at the front of the magnet for patient access and above the magnet for servicing
4. Fluorescent lighting must not be used in the Magnet Room
5. Lighting must be adjusted using a discrete switch or a variable DC lighting controller
6. SCR dimmers or rheostats must not be used
7. LED lighting may be used if the power source is located outside the Magnet Room RF Shield
8. Battery chargers (e.g., used for emergency lighting) must be located outside the Magnet Room
9. Short filament length bulbs are recommended
10. Linear lamps are not recommended due to the high burnout rate

### 9.3 Grounding Requirements

1. The Penetration Panel Pen Wall and Secondary Pen Wall (SPW) are connected to the RF Common Ground Stud with the GE Healthcare supplied ground cable (refer to [Illustration 3-14](#) [Chapter 4, PEN and SPW Wall Opening Requirements](#) )
2. RF Power Filters over 30 volts must be located within 40 in. (1016 mm) of the RF Common Ground Stud
3. RF Power Filters of 30 volts or less may be located anywhere on the RF Shield
4. All electrical devices (e.g., outlets, light fixtures, etc.) must have a ground wire from device power source and be grounded to the RF Shield at the RF Common Ground Stud
5. Resistance between any two grounded devices must not exceed 0.1 ohm to ensure equal potential ground system within the Magnet Room (e.g., MGD to PDU)
6. Do not ground non-MR equipment to the MR ground system
7. All metallic pipes (including water, medical gas, sprinklers, etc.) entering the RF Shield, excluding the Cryogenic Vent and floor drains, must be located within 56 inches (1422 mm) of the RF Common Ground Stud
8. The illustration below shows a typical ground layout

Illustration 3-14: Typical Magnet Room Grounding

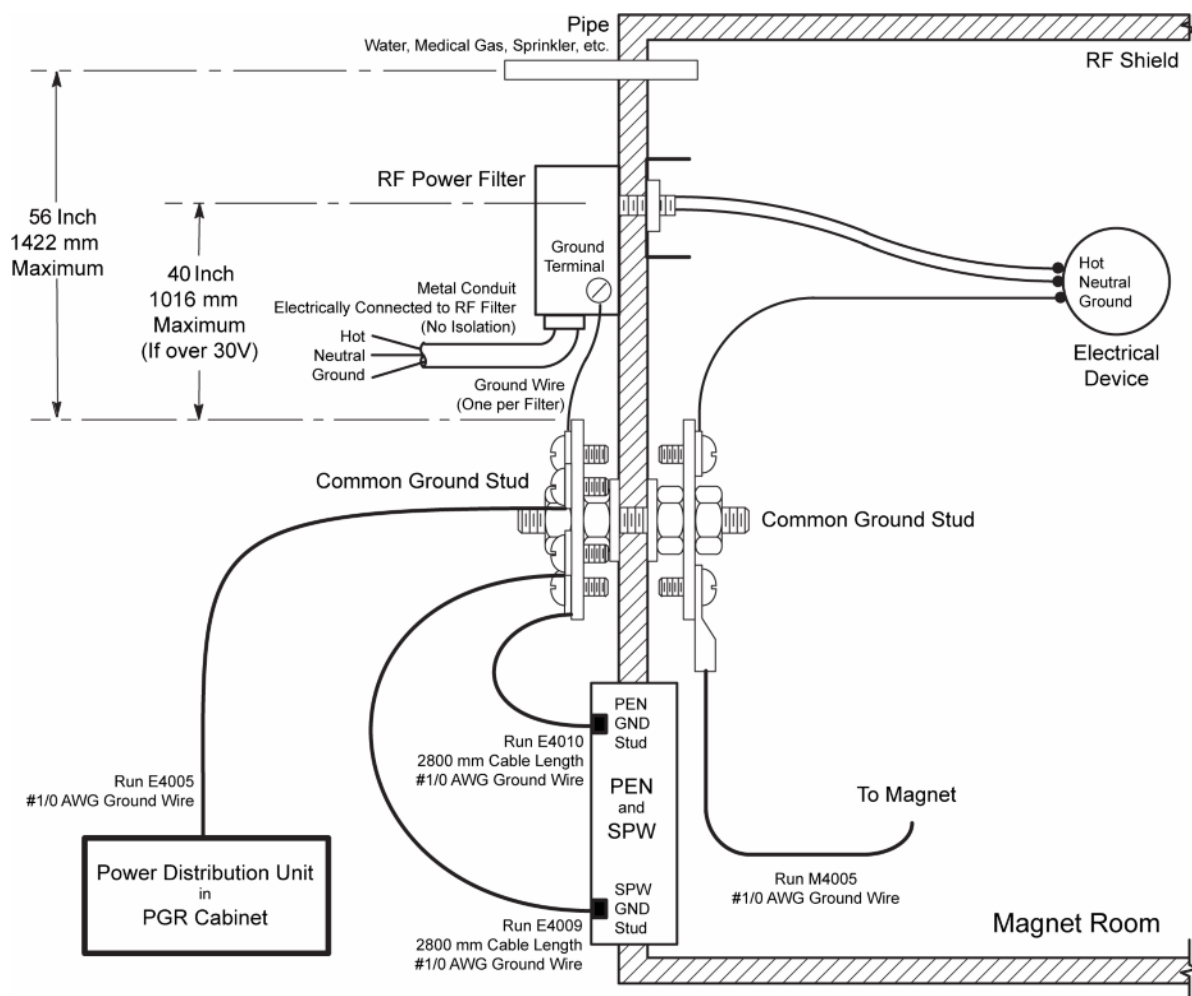


**NOTE:** See MR System Interconnects Specifications for usable cable lengths.

# Chapter 4 Equipment Room

## 1 Equipment Room Overview

The illustration below shows a typical Equipment room layout.

**Illustration 4-1: Typical Equipment Room Layout**

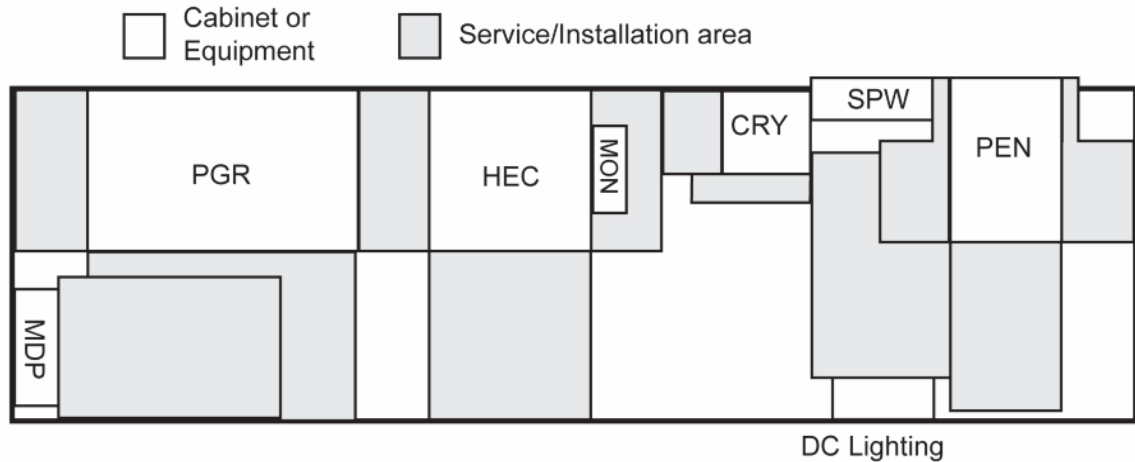

**Notes:**


- Layout shown is typical for an approximate 6 ft (1.8) x 20 ft (6 m) Equipment room
- Drawing is to scale, but see individual component descriptions and room requirements for dimensions
- MDP and Lighting converter service access area long dimension shown as 48 in (122 cm)
- Service/Installation areas may overlap as necessary
- Magnet Monitor (MON) may be mounted on either side of the HEC
- Optional equipment is not shown

## 2 Main Disconnect Panel (MDP) Specifications

The Main Disconnect Panel (MDP) is provided with the MR system.

1. Weight: 130 lbs (59 kg)
2. Magnetic Field Limit: 50 Gauss (5 mT)


**Illustration 4-2: Main Disconnect Panel**

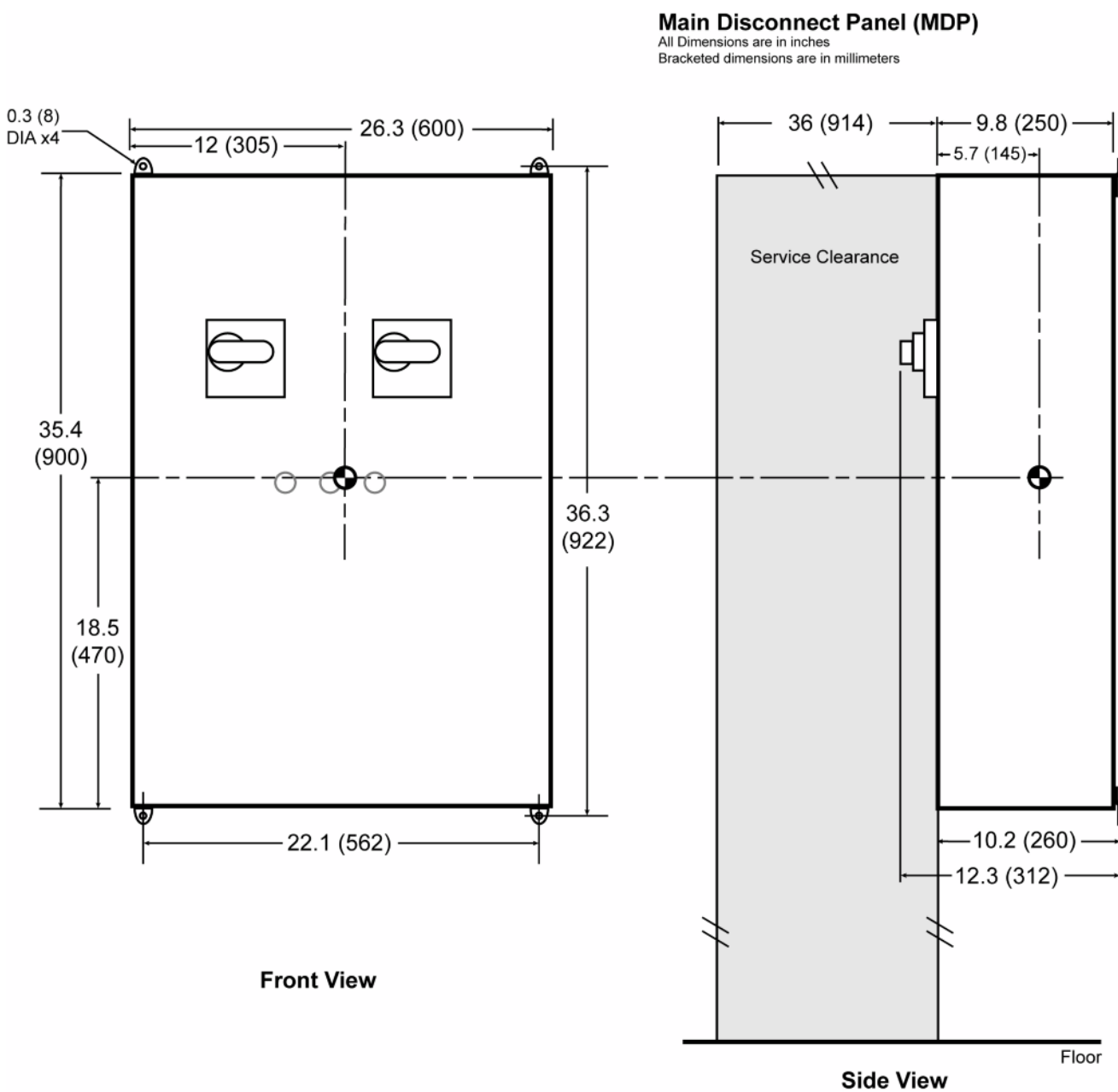


### 3 Power, Gradient, RF Cabinet (PGR) Specifications

1. Weight: 3144 lbs (1426 kg)
2. Magnetic Field Limit: 50 Gauss (5 mT)

Illustration 4-3: Power, Gradient, RF (PGR) Cabinet

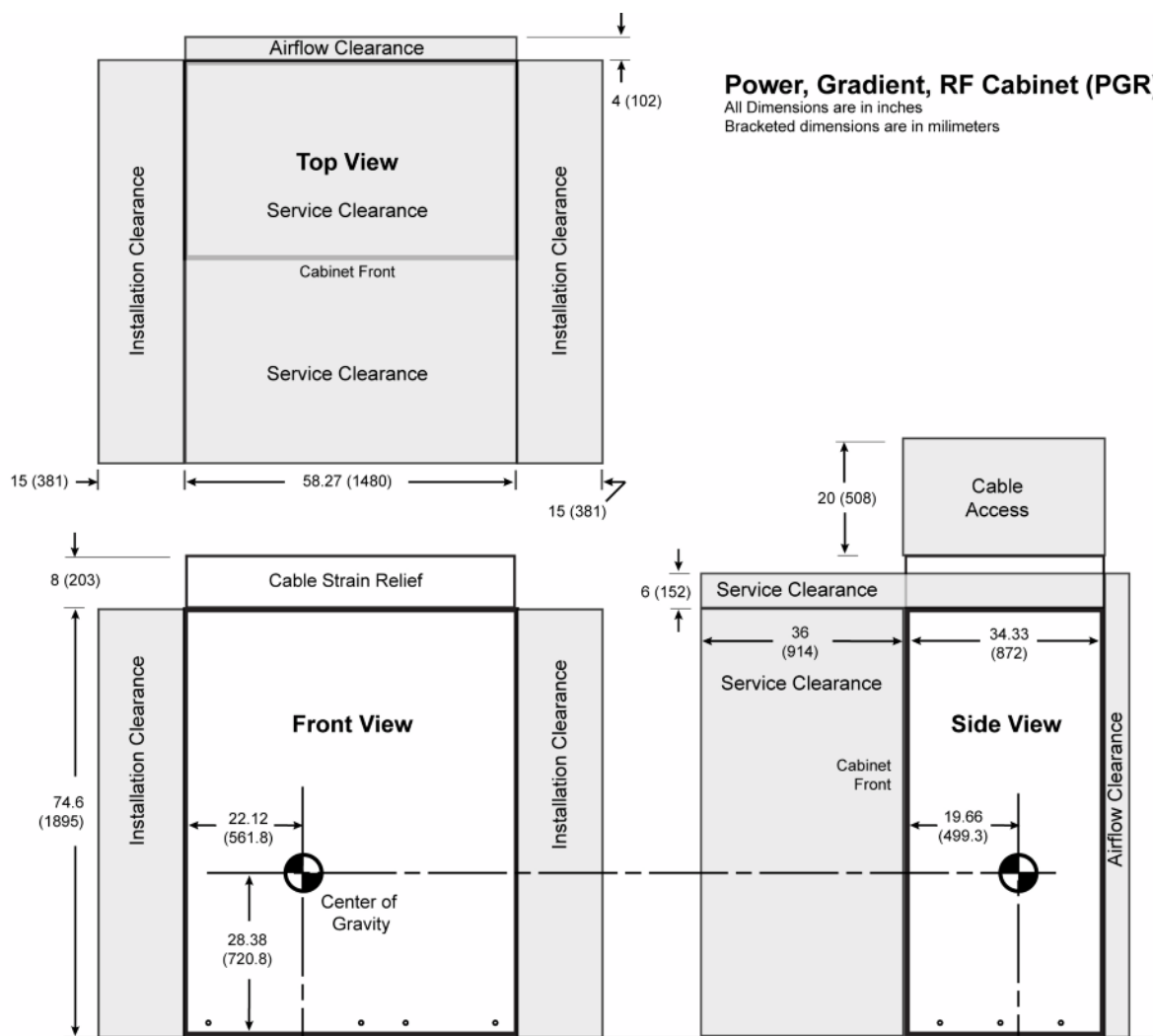

#### Mounting Hole Detail

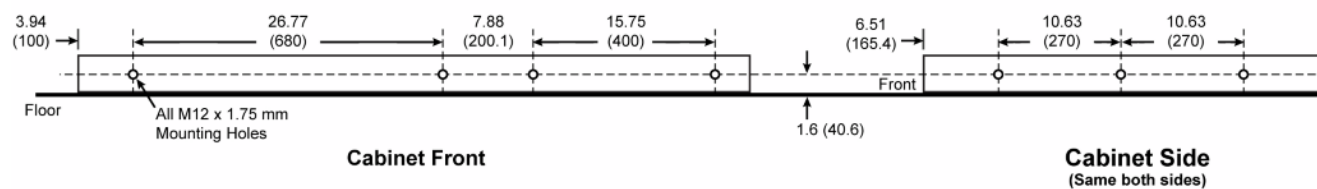


## 4 Heat Exchanger Cabinet (HEC) Specifications

Facility coolant is supplied to the Heat Exchanger Cabinet (HEC). The HEC provides closed loop chilled water circuits for the Power, Gradient, RF (PGR) Cabinet in the equipment room and the Gradient Coil inside the Magnet Enclosure. The HEC also routes facility chilled coolant to the Cryocooler Compressor (CRY).

1. Weight (approximate):
  - a. Dry (shipping – no fluid): 950 lbs (431 kg)
  - b. Wet (with cooling fluid): 1350 lbs (612 kg)
2. Magnetic Field Limit: 50 Gauss (5 mT)



Illustration 4-4: Heat Exchanger Cabinet (HEC)

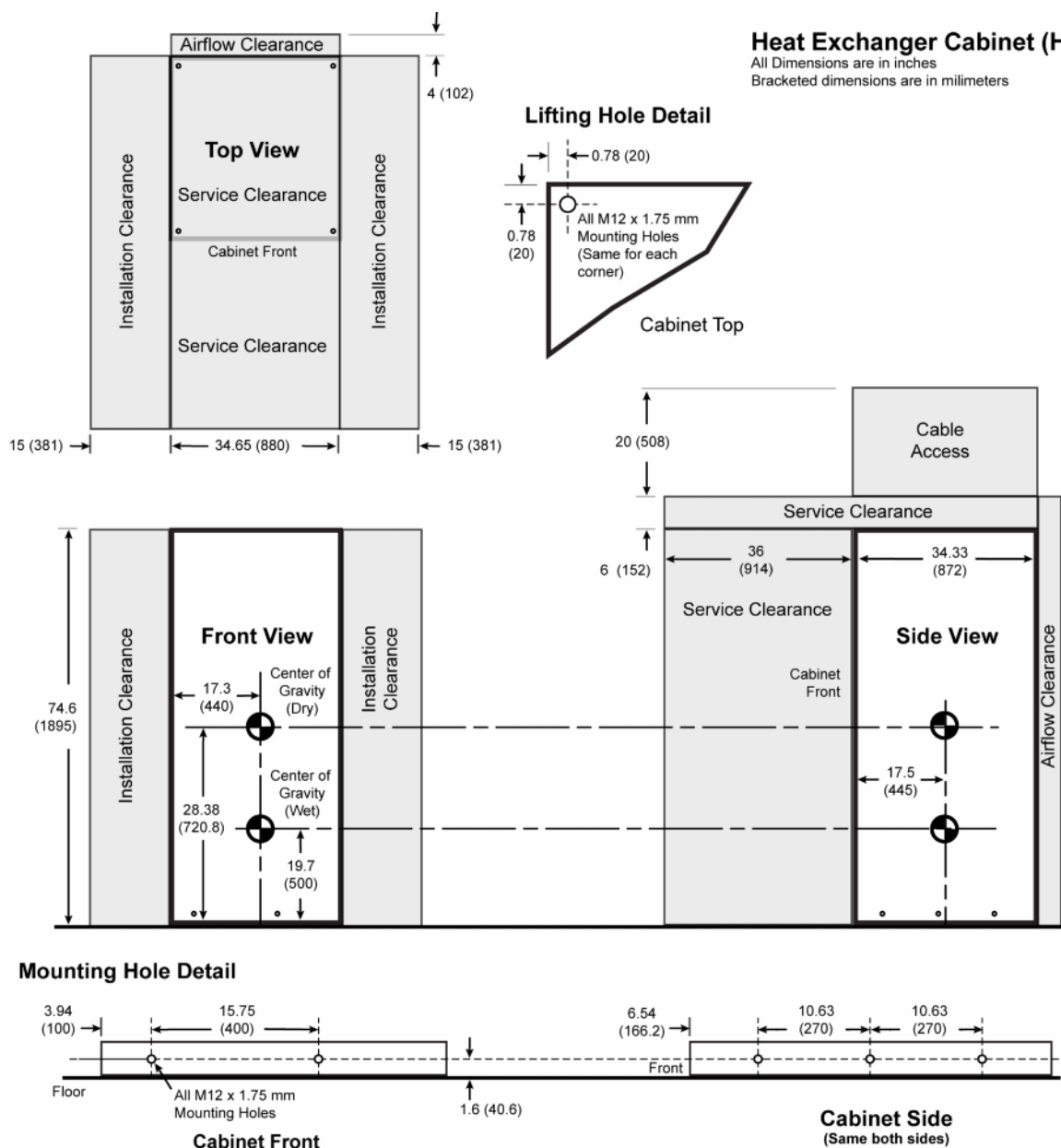


## 5 Penetration Panel Cabinet (PEN) Specifications

The Penetration Cabinet Penetration Panel provides interconnects from the PEN Cabinet through the Magnet Room RF Shield.

1. Weight: 639 lbs (290 kg)
2. PEN Panel Magnetic Field Limit: 200 Gauss (20 mT) for the entire PEN Panel (i.e., the blower box must be outside the 200 Gauss line)  
PEN Cabinet Magnetic Field Limit: 50 Gauss (5 mT)
3. The Penetration Panel Blower Box in the Magnet room requires airflow of 120 CFM with a maximum pressure drop of 0.5 inches of water (gage)
4. The PEN cabinet must be positioned directly in front of the PEN Panel. See [PEN and SPW Wall Opening Requirements](#) for PEN panel mounting and location requirements

Illustration 4-5: Penetration Panel Cabinet (PEN)

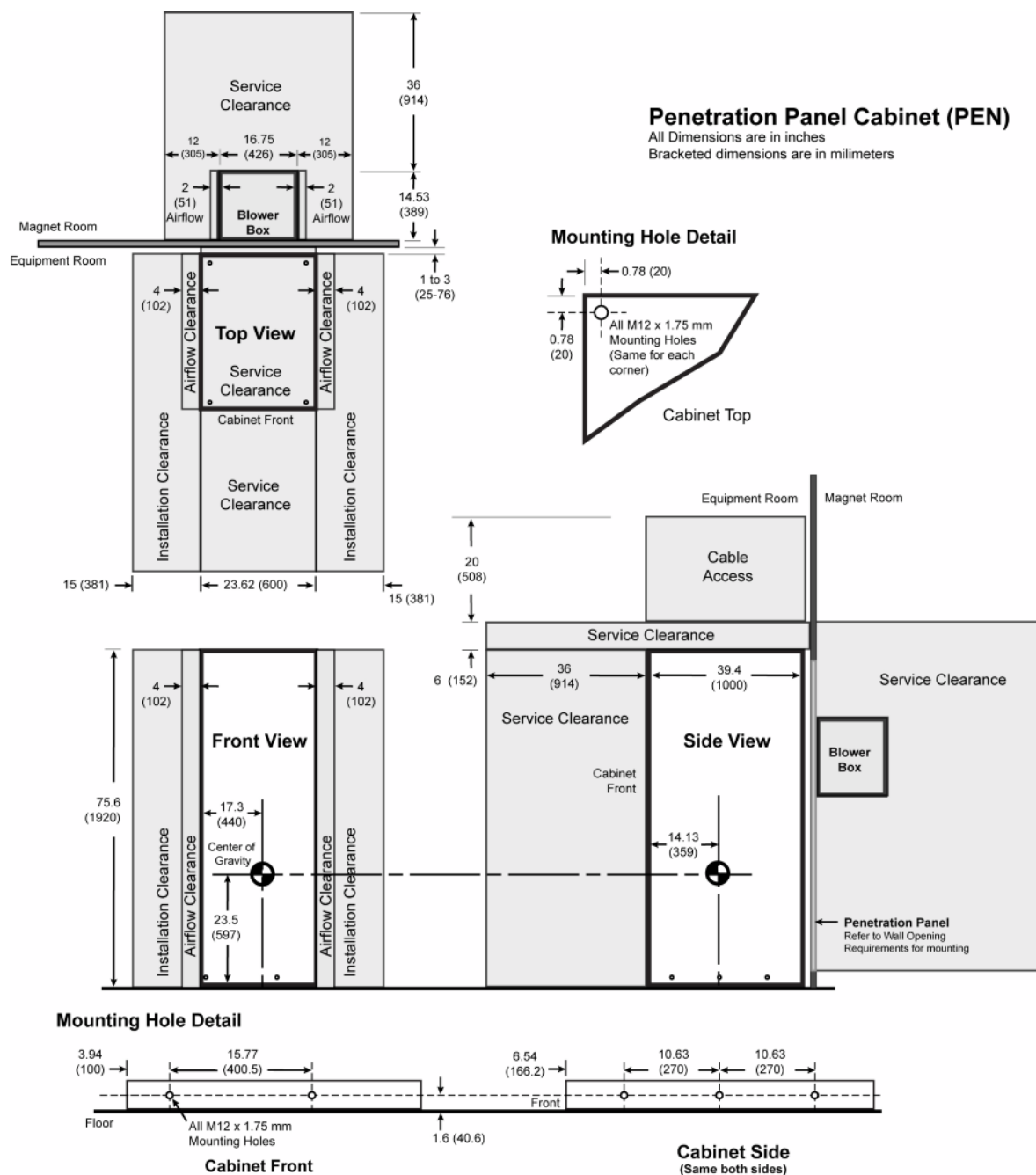


## 6 Secondary Penetration Wall (SPW) Specifications

The Secondary Penetration Wall (SPW) provides interconnects from the Equipment Room through the Magnet Room RF Shield.

1. Maximum Magnetic Field: 200 Gauss (20mT)
2. See [PEN and SPW Wall Opening Requirements](#) mounting and location requirements


### Illustration 4-6: Secondary Penetration Wall (SPW)

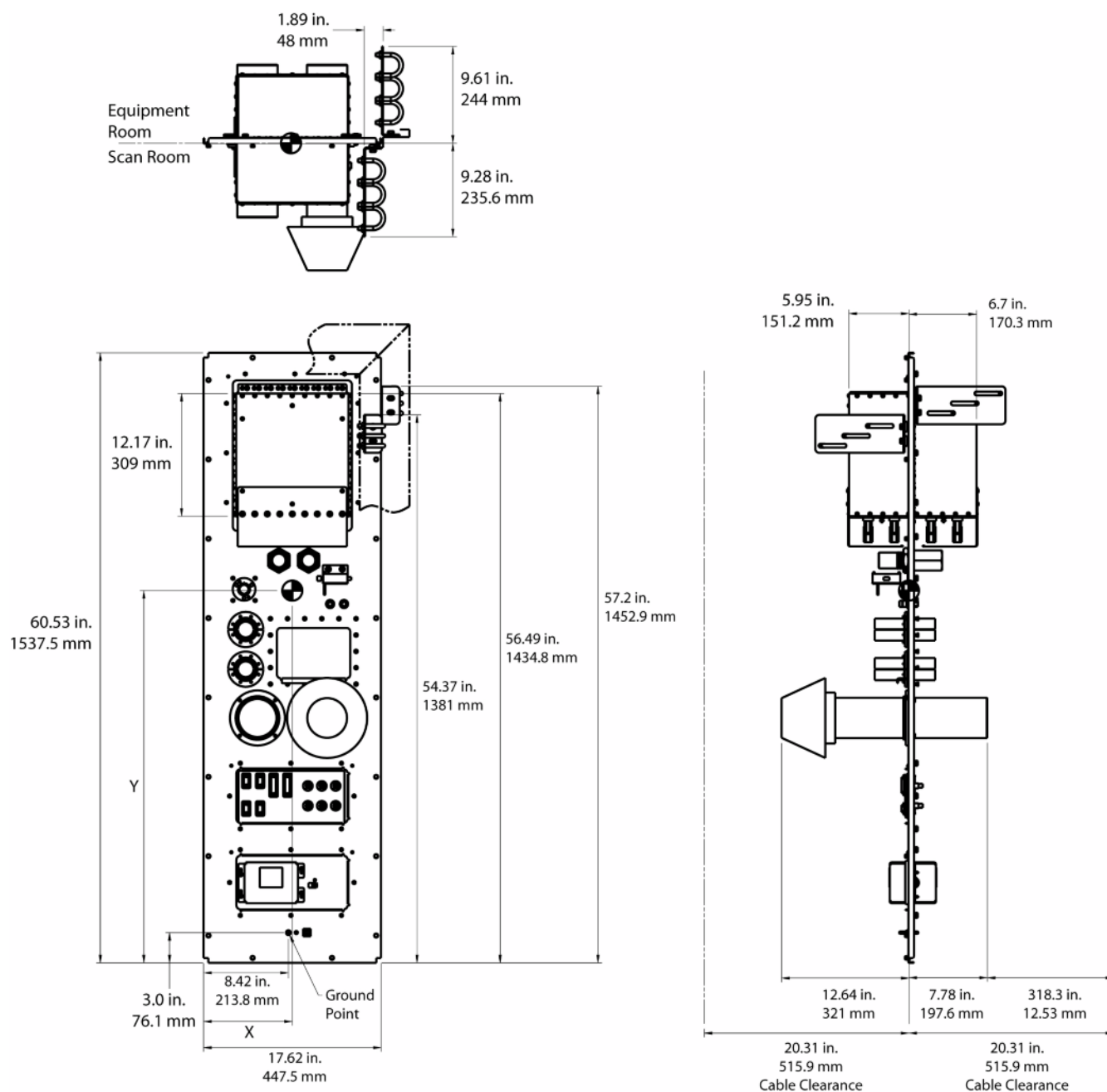

## 7 PEN and SPW Wall Opening Requirements

1. The Equipment Room and the Magnet Room must share at least one common wall to mount the PEN panel and SPW
2. The penetration panel opening requirements are shown in [Illustration 4-7](#)

**Illustration 4-7: PEN and SPW Panel Openings**

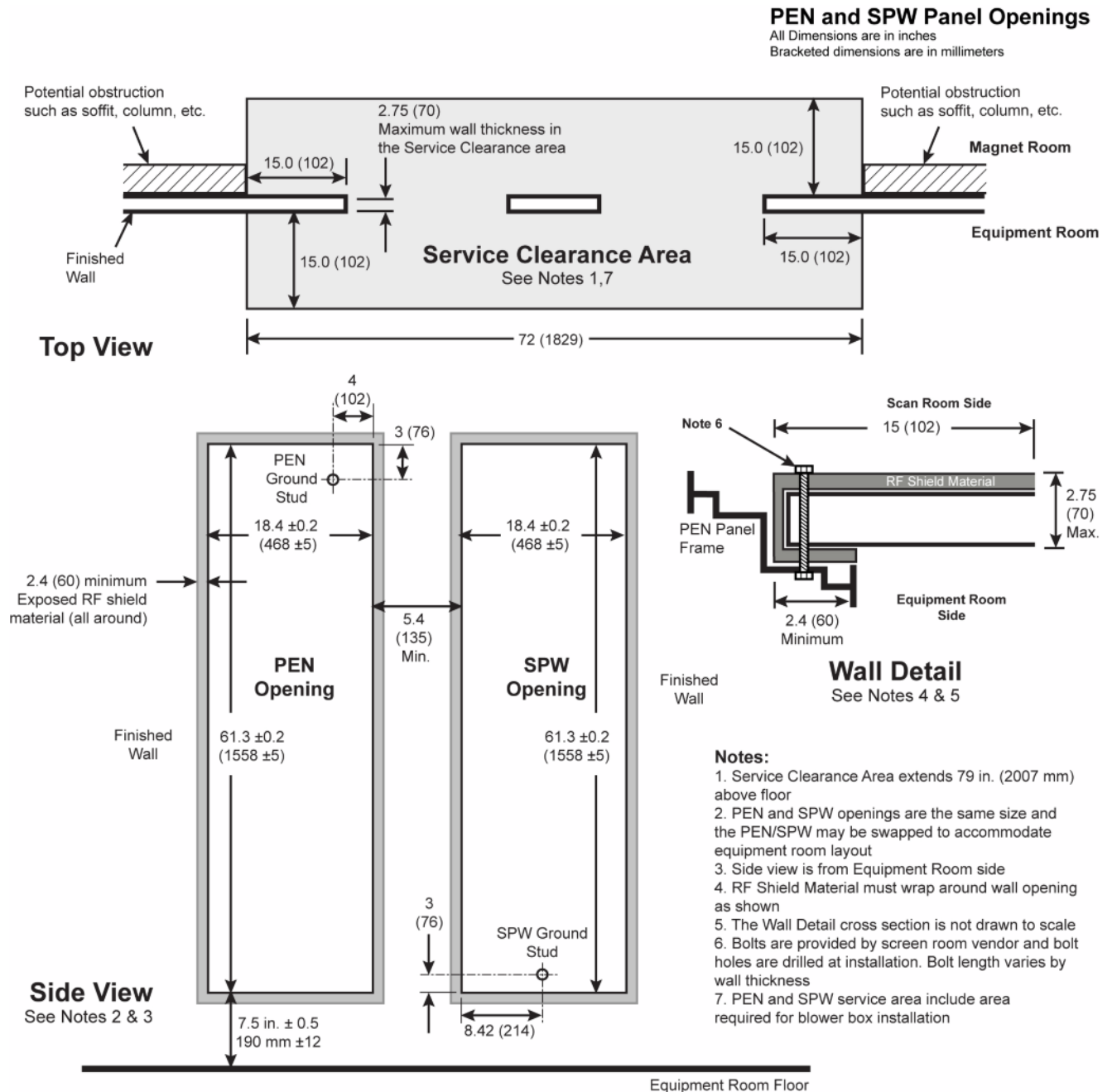


## 8 Cryocooler Compressor (CRY) Specifications

Water cooling for the Cryocooler Compressor (CRY) is provided from the Heat Exchanger Cabinet (HEC) or facility supplied emergency backup water supply.

1. Weight: 264 lbs (120 kg)
2. Magnetic Field Limit: 100 Gauss (10 mT)

**Illustration 4-8: Cryocooler Compressor F-50**

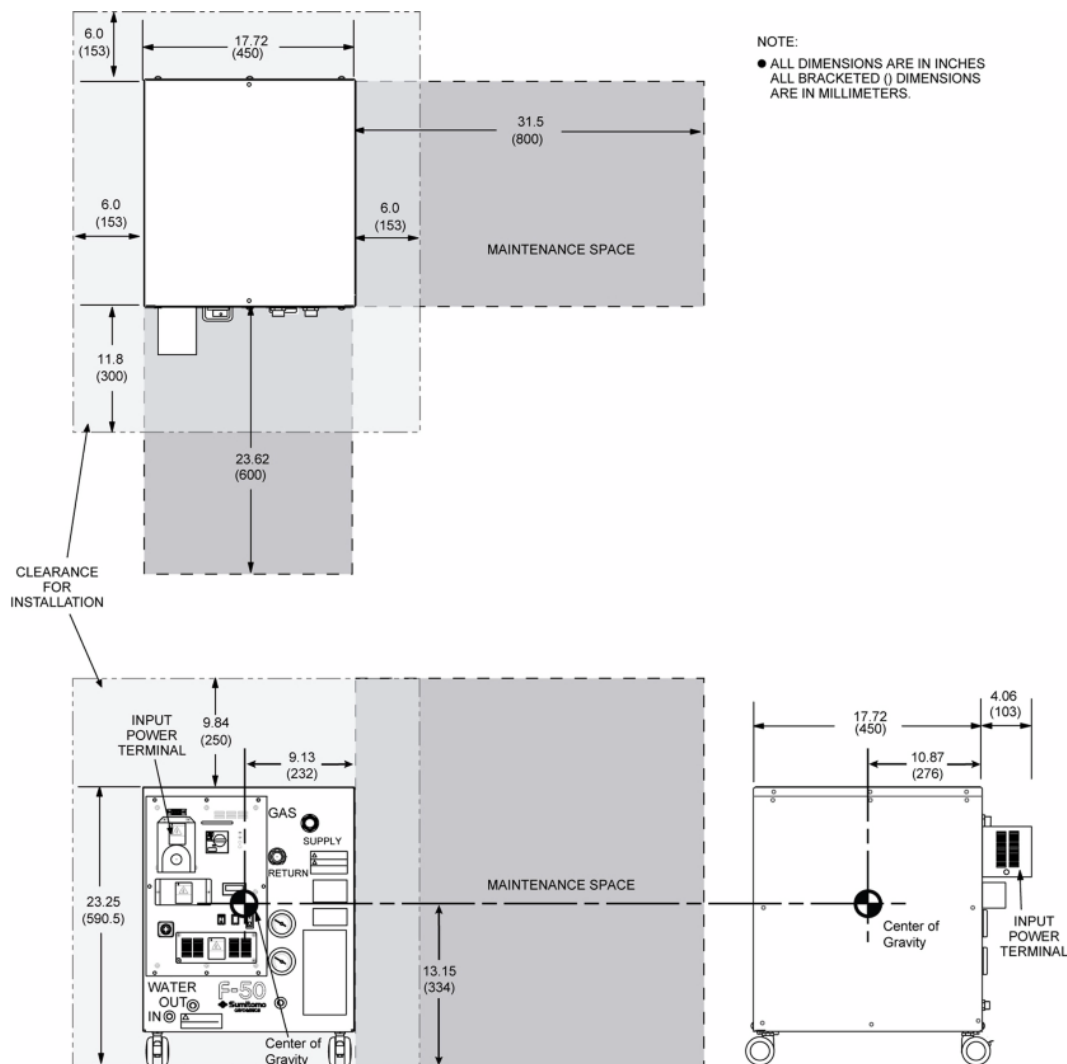

## 9 Magnet Monitor (MON) Requirements and Specifications


### 9.1 Requirements

1. Customer must supply T100 network connection with RJ45 connector to the Magnet Monitor (MON)
2. The network connection must not be routed through the Ethernet switch in the Global Operator Cabinet (GOC)
3. The Magnet Monitor requires a 110/220 VAC, 50/60 Hz, 2.0 A facility supplied outlet. Power at the outlet must be continuously available

### 9.2 Specifications

1. Mounting location: On either side of the Heat Exchanger Cabinet (HEC)
2. Weight: 10 lbs (4.5 kg)
3. Maximum gauss limit: 200 gauss (20 mT)
4. Power cord length: 72 in. (1829 mm)

Illustration 4-9: Magnet Monitor (MON)

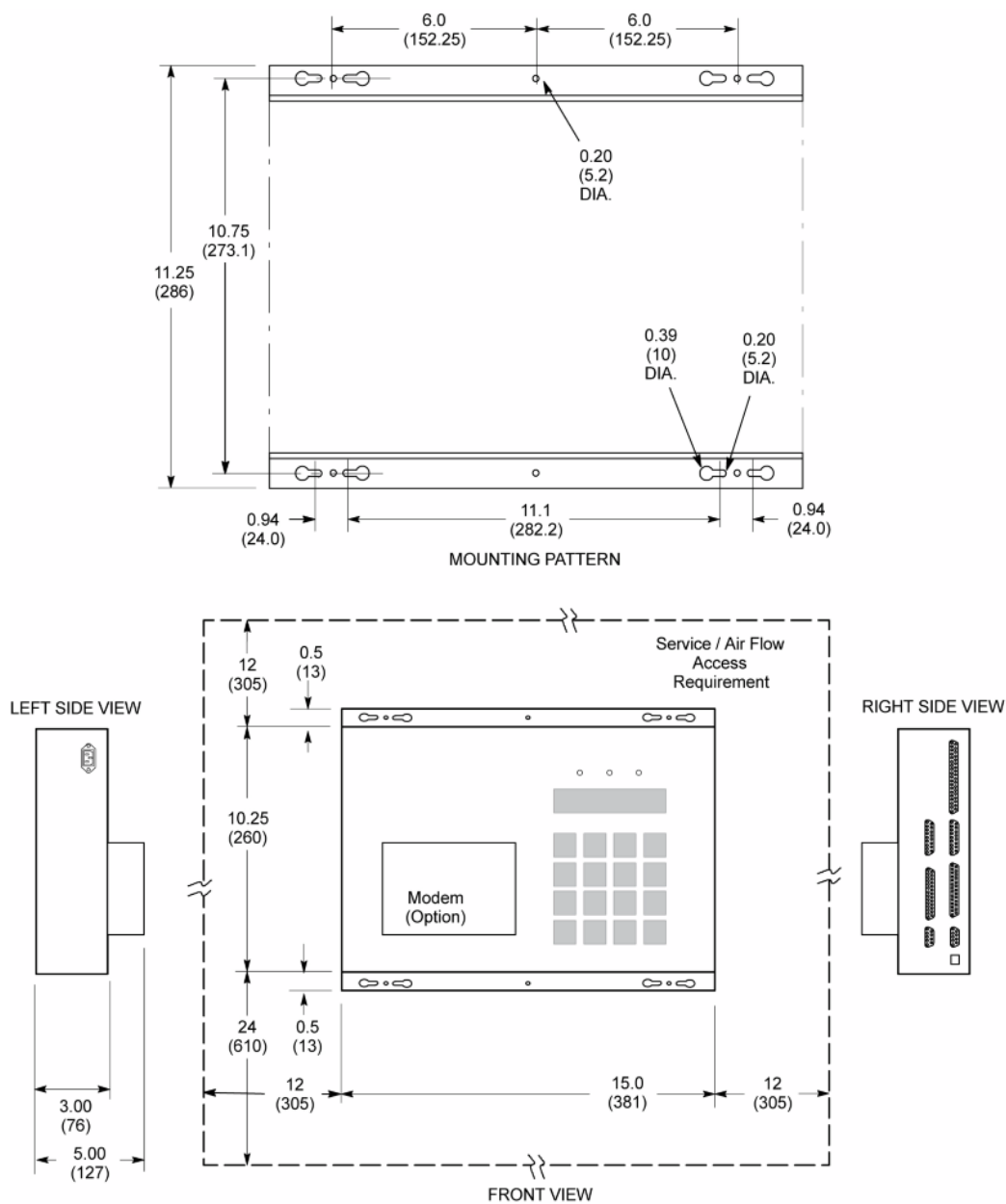


All dimensions are in inches, all bracketed ( ) dimensions are in millimeters.



## 10 Brainwave Lite (BW) Specifications

Illustration 4-10: Brainwave Lite Cabinet (BW)

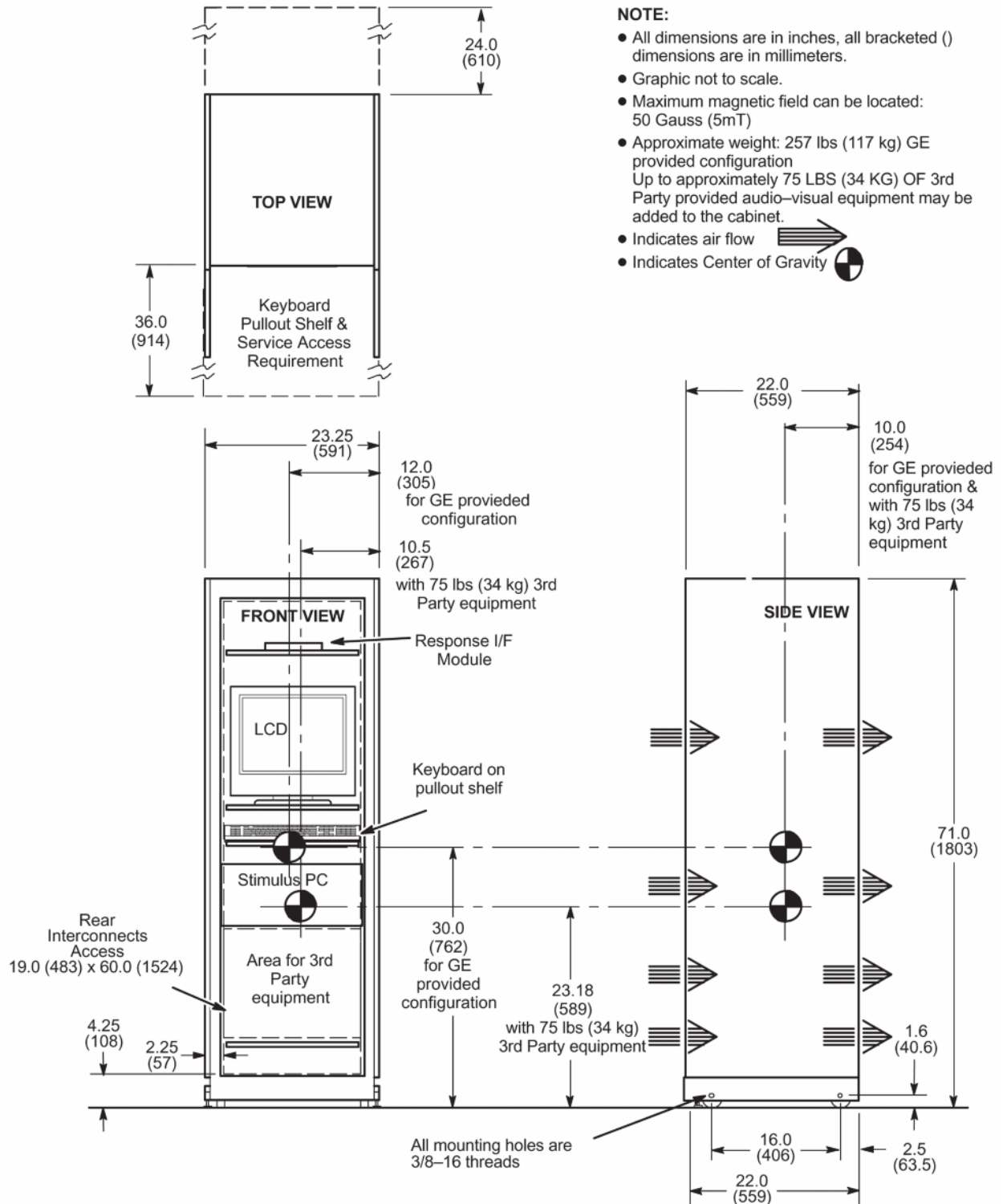

## 11 Magnetic Resonance Elastography (MRE) Specifications


### 11.1 Requirements

1. Customer must work with the RF shield vendor to provide a waveguide for the 1 in. (25.4 mm)
2. MRE Resonant Acoustic Driver location is limited to the length of the 1 in. (25.4 mm) tube (see the usable cable lengths in [Chapter 6, Section 1.4, Magnetic Resonance Elastography \(MRE\) Option](#))
3. The MRE Resonant Acoustic Driver module requires a 110/220 VAC, 50/60 Hz, 2.0 A facility supplied outlet. Power at the outlet must be continuously available

### 11.2 Specifications

1. Weight: 53.4 lbs (24.22 kg)
2. Gauss Limit: 50 gauss
3. Power Cord Length:
  - 60 Hz: 240 in. (6096 mm)
  - 50 Hz: 300 in. (7620 mm)

**Illustration 4-11: Magnetic Resonance Elastography (MRE) Resonant Acoustic Driver**

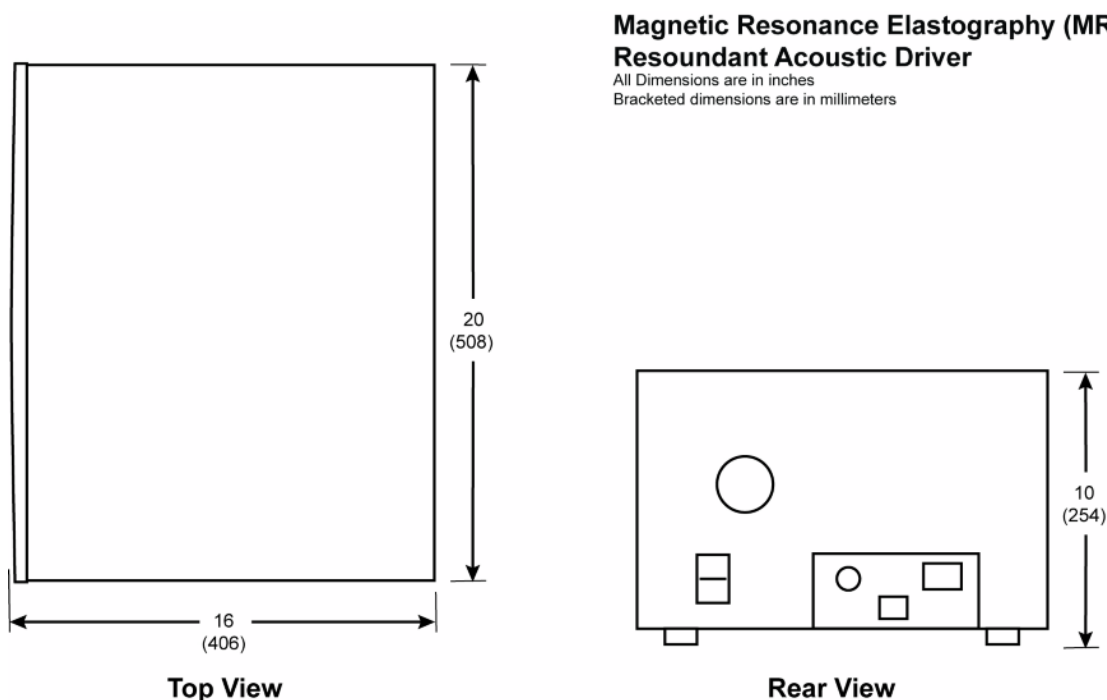


## Chapter 5 Control Room

### 1 Operator Workspace Equipment Specifications

The operator seated at the Operator Workspace must have an unobstructed view of the patient on the Patient Table docked to the Magnet.

#### 1.1 Operator Workspace (OW)

Illustration 5-1: Operator Workspace (OW)

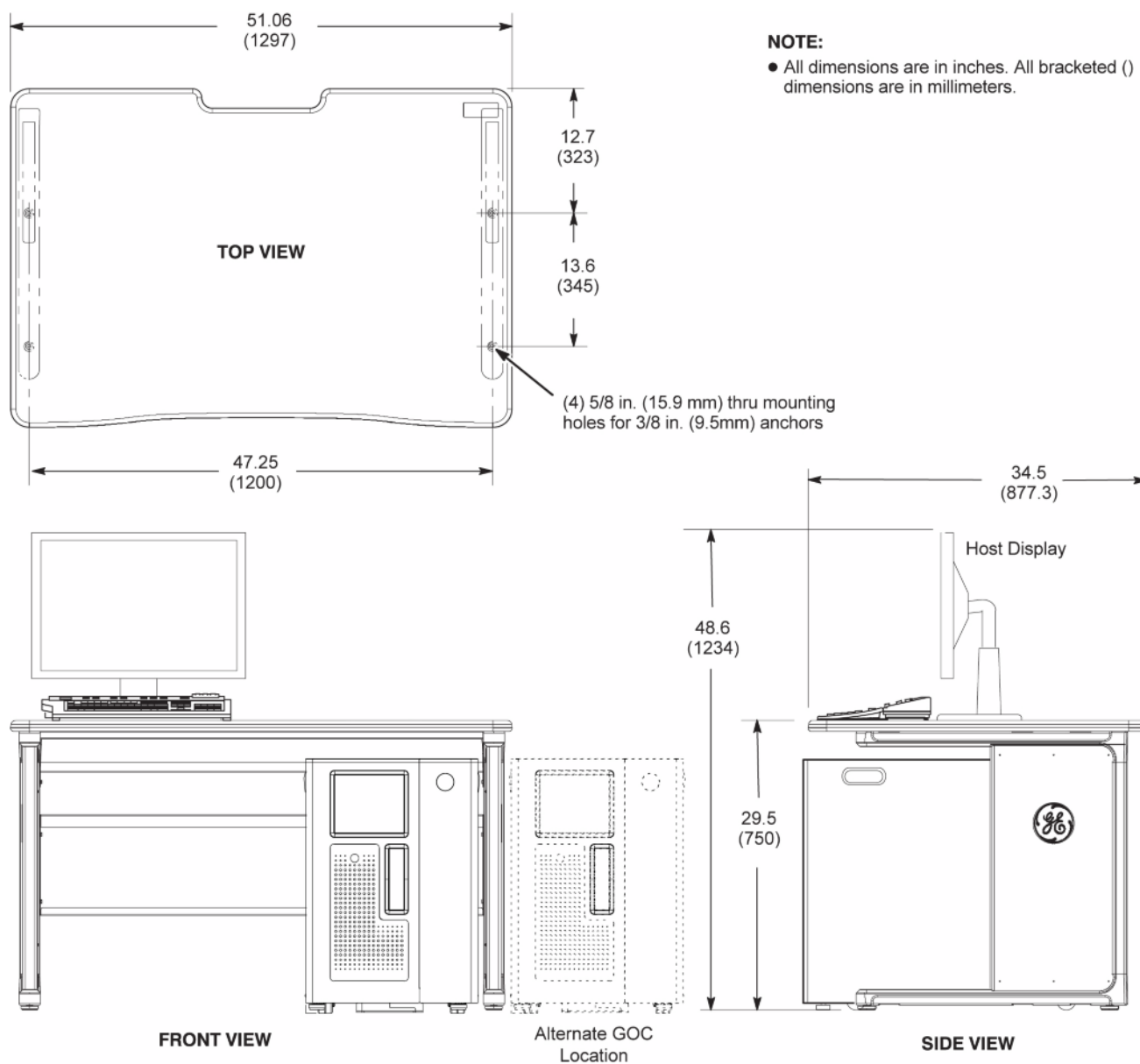


#### 1.2 Global Operator Cabinet (GOC)

1. Weight: 141.75 lbs (64.3 kg)

## 2. Magnetic Field Limit: 50 Gauss (5 mT)


Illustration 5-2: Global Operator Cabinet (GOC)

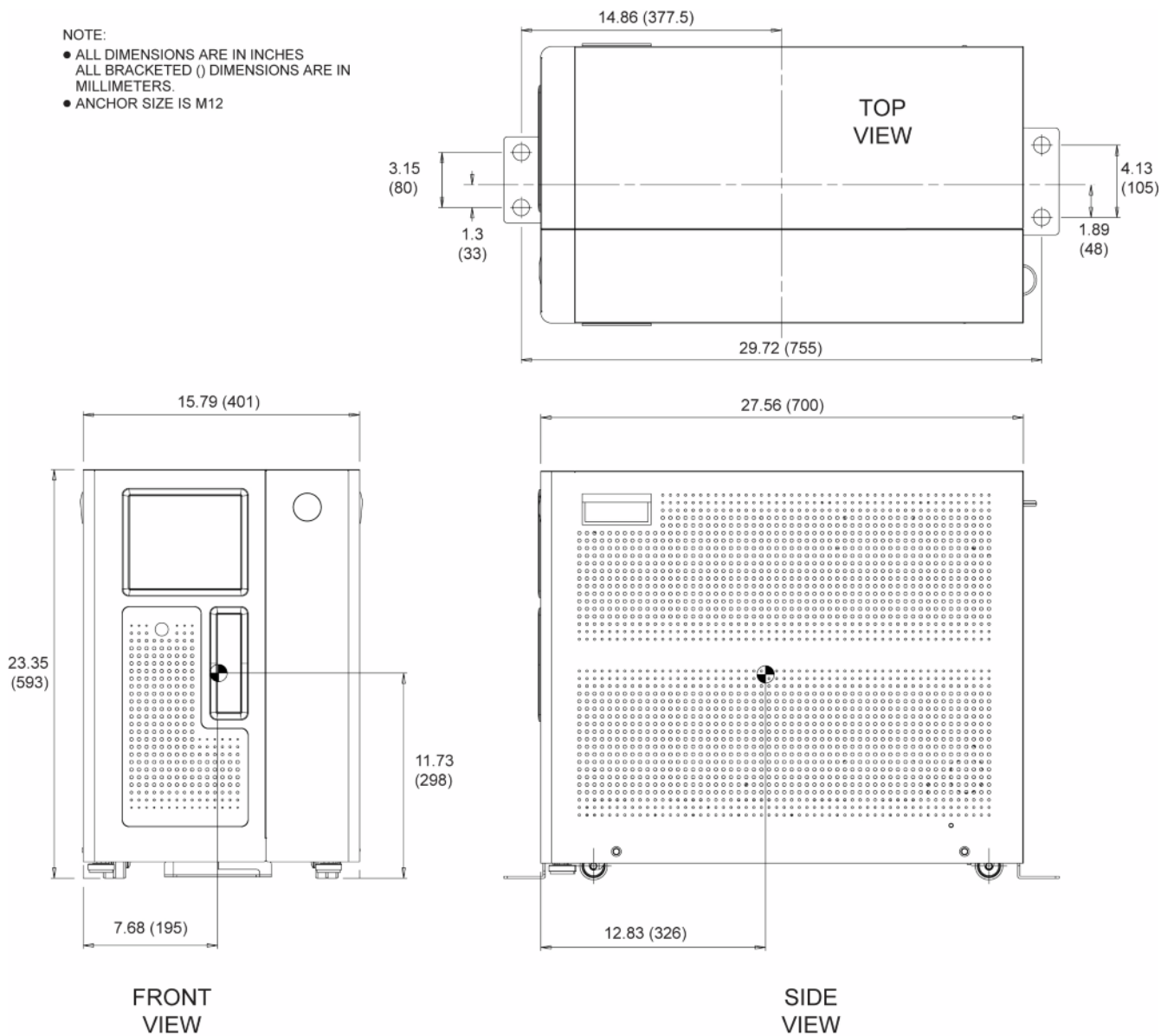

## 1.3 Host Display


1. Weight 26 lbs (11.8 kg)
2. Magnetic Field Limit: 50 Gauss (5 mT)

Illustration 5-3: Host Display

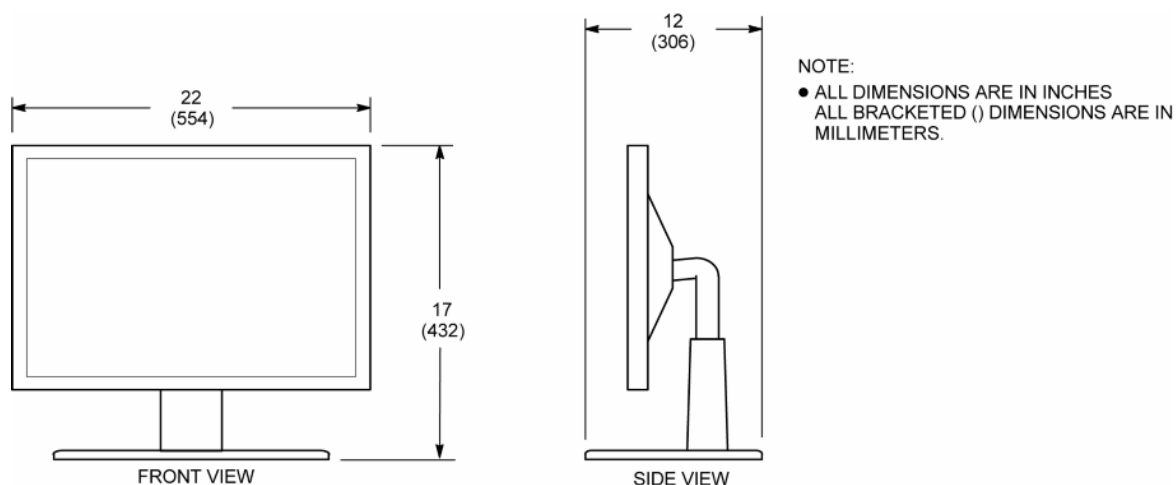

## 1.4 Host Keyboard

Illustration 5-4: Host Keyboard

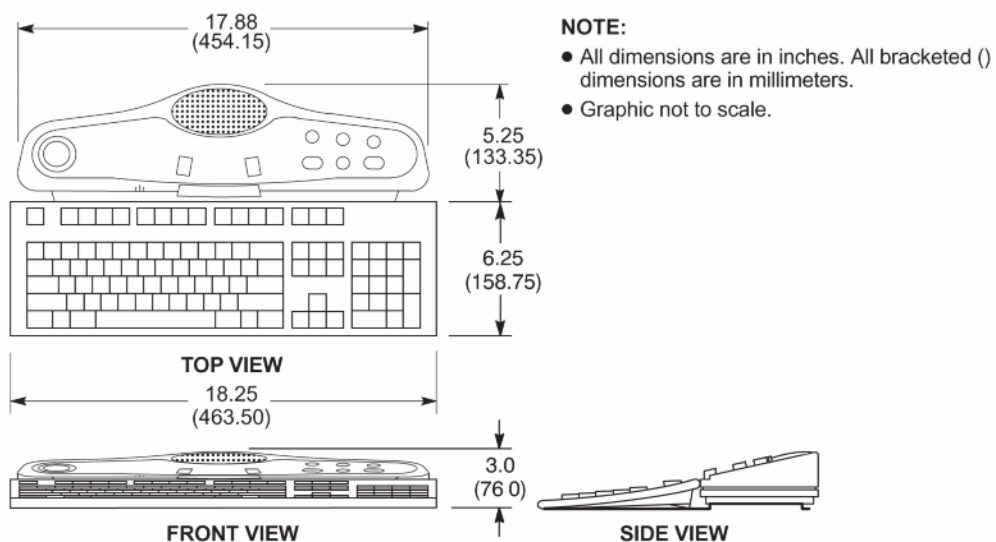


## 1.5 Pneumatic Patient Alert

The Pneumatic Patient Alert system allows the patient to contact the operator. The Control Box audible and visual alarm will be activated by the patient squeeze bulb which is located on the Magnet Enclosure and connected by pneumatic tubing through the Penetration Panel to the Control Box.

1. Weight 0.5 lbs (0.2 kg)
2. Magnetic Field Limit: 50 Gauss (5 mT)
3. The Control Box must be mounted within sight of operator and within 5 ft. (1.5 m) of an electrical outlet

Illustration 5-5: Pneumatic Patient Alert (PA)

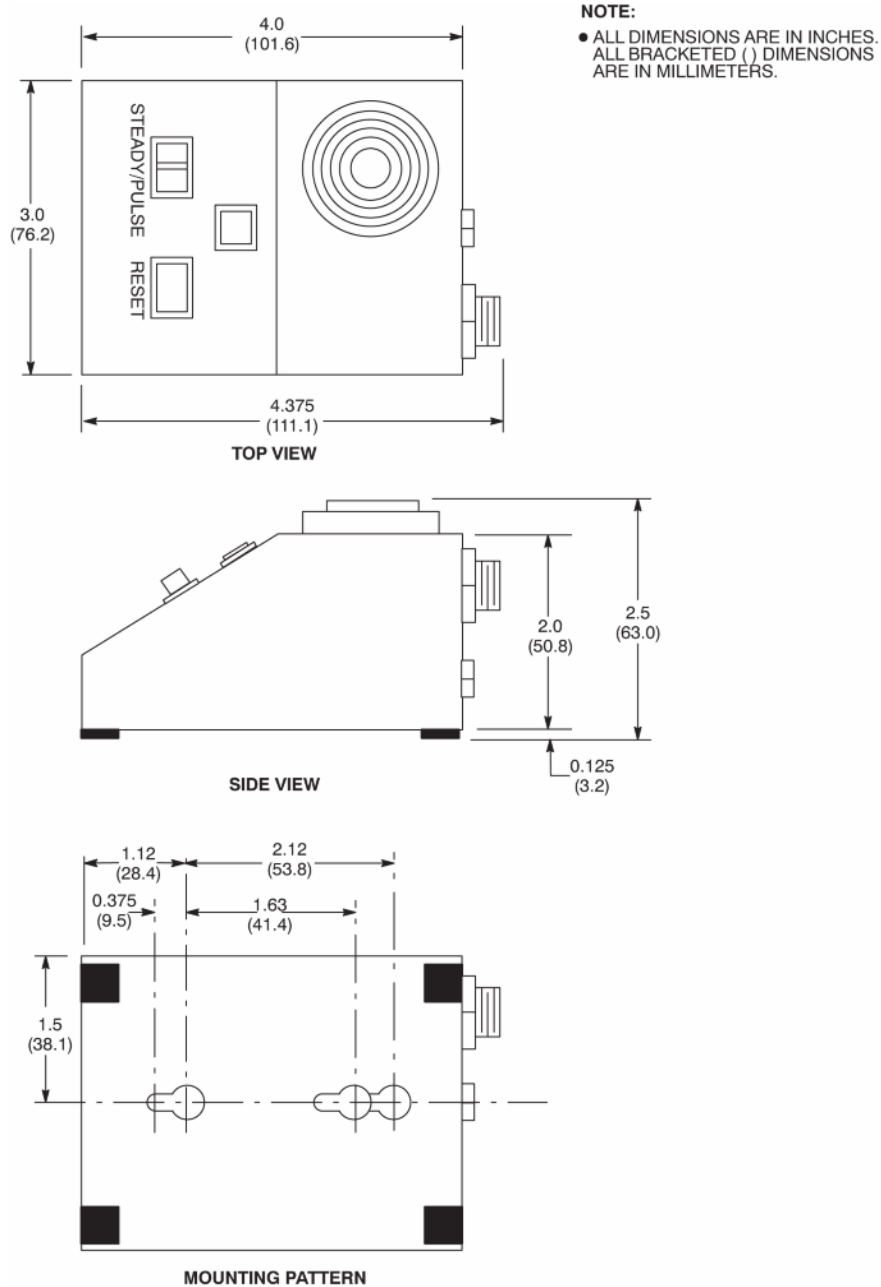

## 2 Oxygen Monitor (OXY) Specifications

The Oxygen Monitor system consists of the Oxygen Monitor, the Remote Oxygen Sensor Module, and interconnects through the Secondary Penetration Wall (SPW). The Oxygen Monitor alarm is located near the Operator Workspace is activated by the Remote Oxygen Sensor Module in the Magnet Room.

### 2.1 Oxygen Monitor (OXY) Requirements

The Oxygen Monitor requires facility supplied power


**Table 5-1: Oxygen Monitor Facility Power Requirements**

| Parameter           | Requirements       |          |
|---------------------|--------------------|----------|
| Voltage / Frequency | 100-120 VAC        | 50/60 Hz |
|                     | 200-220 VAC        | 50/60 Hz |
| Phase               | 1                  |          |
| Maximum Amps        | 3.0                |          |
| Connection type     | Hard wired in unit |          |

Illustration 5-6: Oxygen Monitor (OXY)

**NOTE:**

- All dimensions are in inches. All bracketed dimensions are in millimeters.
- Graphic not to scale.
- Maximum magnetic field can be located: 100 Gauss (10mT)
- Approximate Weight: 9 lbs (4.1 kg)

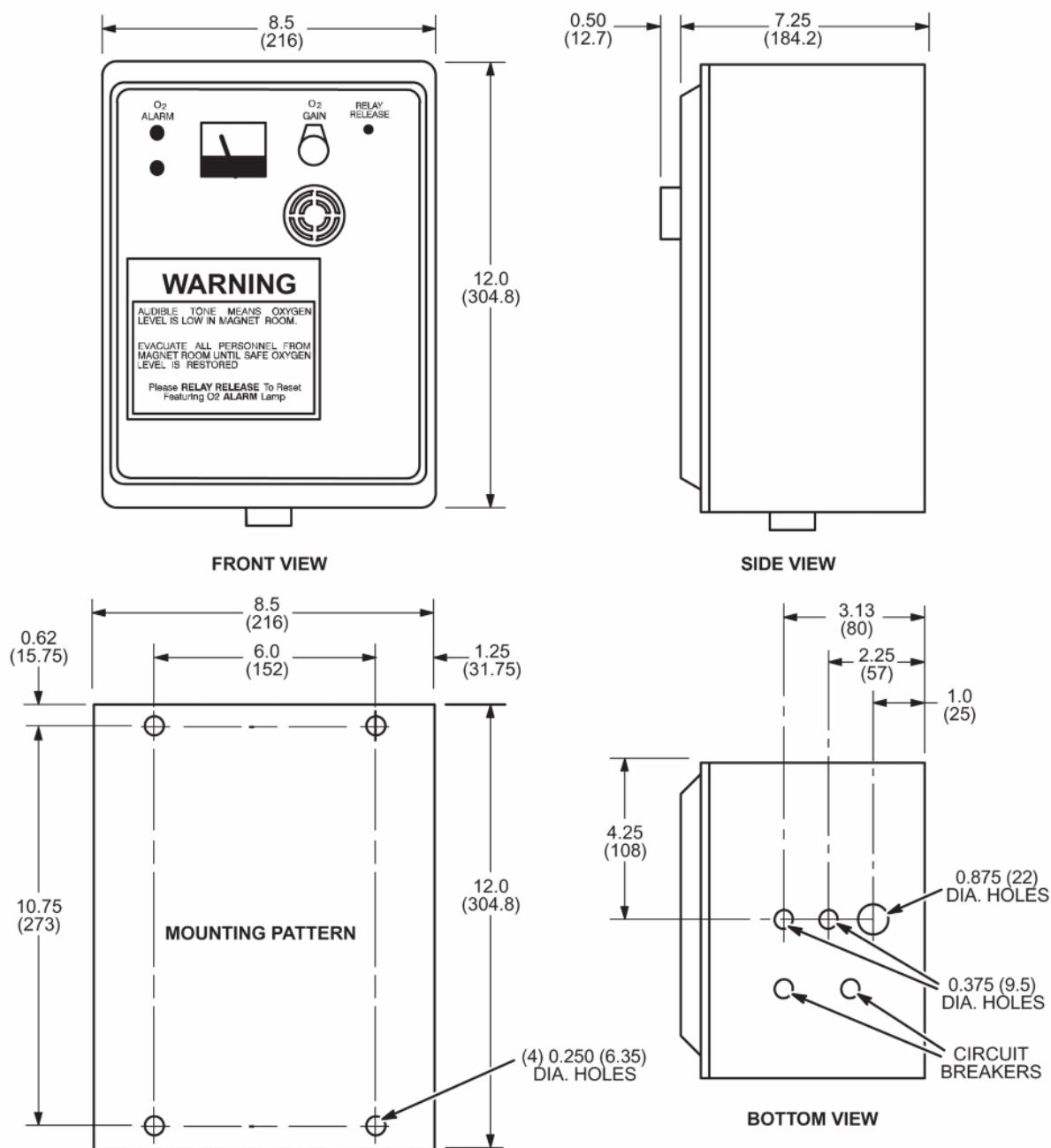




## **2.2 Oxygen Monitor Sensor**

The Remote Oxygen Sensor Module must be mounted approximately 60 in. (1524 mm) above the Magnet Room floor near the front of the Magnet Enclosure.

### Illustration 5-7: Oxygen Monitor Sensor

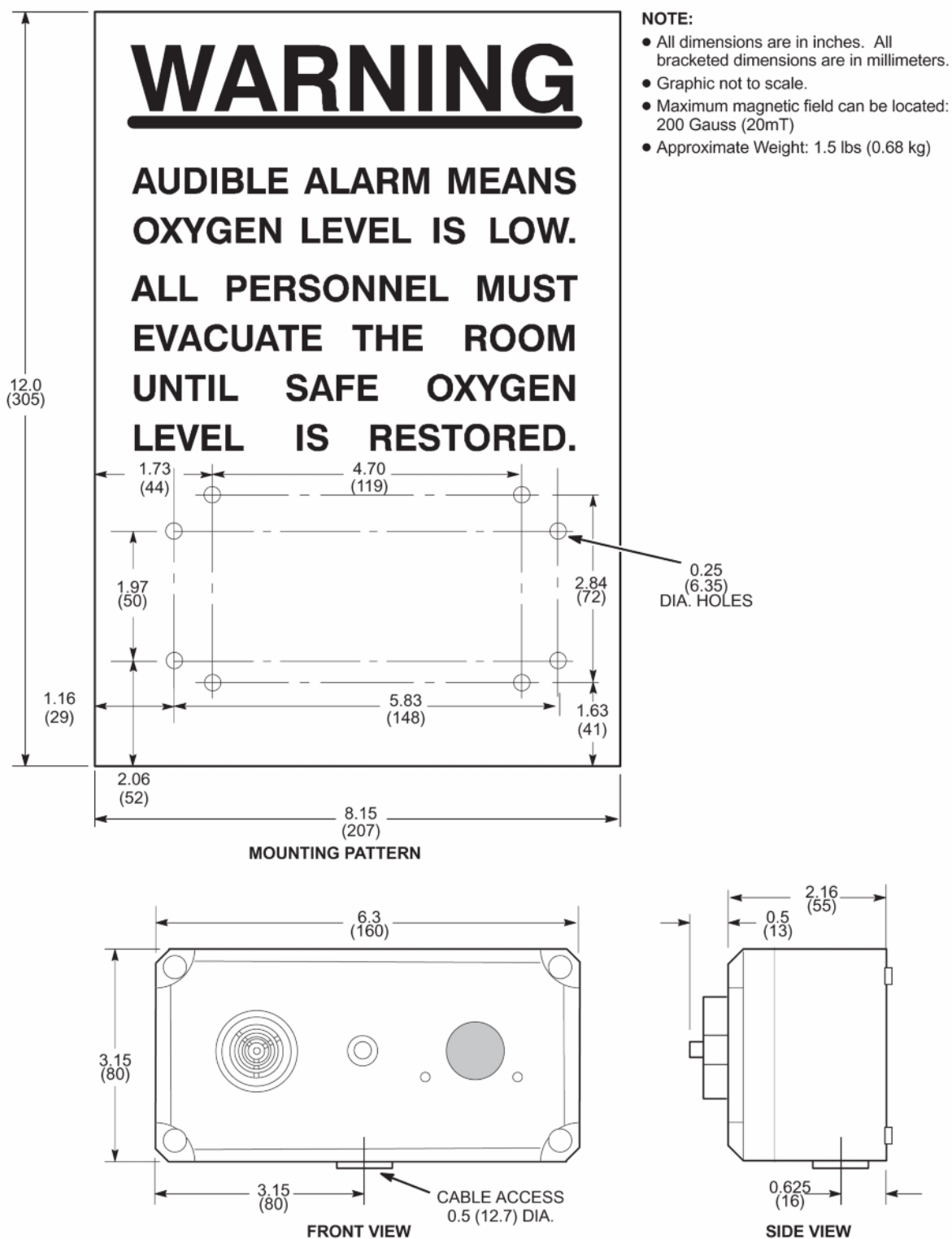

# Chapter 6 MR System Interconnects

## 1 MR System Interconnects Specifications

### 1.1 Component Designator Definitions

GE Healthcare uses Component Designators to identify system components. All subsystem cabinets and other components are referred to by their component designators in the Interconnect Data diagrams and tables.

**Table 6-1: MR System Component Designators**

| Component Designator | Description                                                           |
|----------------------|-----------------------------------------------------------------------|
| CRY                  | Cryocooler Compressor Cabinet                                         |
| E01/E02              | Emergency Off Buttons                                                 |
| MDP                  | Main Disconnect Panel                                                 |
| HEC                  | Heat Exchanger Cabinet                                                |
| MAG                  | Magnet and Enclosure (all magnet enclosure components in Magnet Room) |
| Modem                | Modem for Magnet Monitor                                              |
| MON                  | Magnet Monitor                                                        |
| MRU                  | Magnet Rundown Unit                                                   |
| OW                   | Operator Workstation                                                  |
| PA1                  | Pneumatic Patient Alert Control Box                                   |
| PDU                  | Power Distribution Unit (PDU) is a module in the PGR cabinet          |
| PED                  | Magnet Rear Pedestal                                                  |
| PEN                  | Penetration Panel Cabinet                                             |
| PGR                  | Power Gradient RF Cabinet                                             |
| PT                   | Patient Transport Table                                               |
| DS                   | Door Switch                                                           |
| SPW                  | Secondary Pen Wall                                                    |

**Table 6-2: MR System Options Component Designators**

| Component Designator | Description                     |
|----------------------|---------------------------------|
| BW                   | Brainwave Lite Cabinet          |
| OXY                  | Oxygen Monitor                  |
| OM2                  | Remote Oxygen Sensor Module     |
| MRE                  | Magnetic Resonance Elastography |

### 1.2 Usable Cable Lengths

**Table 6-3: Usable Cable Lengths**

| Point A        | Point B          | Site Option: Short | Site Option: Long |
|----------------|------------------|--------------------|-------------------|
|                |                  | in. (mm)           |                   |
| Equipment Room |                  |                    |                   |
| CRY            | Magnet Monitor   | 370.1 (9400)       |                   |
| Magnet Monitor | SPW, bottom edge | 661.4 (16800)      |                   |

|                        |                                              |                |                |
|------------------------|----------------------------------------------|----------------|----------------|
| HEC, top panel         | CRY                                          | 354.3 (9000)   |                |
| HEC, top panel         | Customer-Supplier Network                    | 1102.4 (28000) |                |
| PEN cabinet, top panel | GOC, rear panel                              | 141.7 (29000)  |                |
| PEN cabinet, top panel | Magnet Monitor                               | 645.7 (16400)  |                |
| PEN Ground stud        | RF common ground stud                        | 86.6 (2200)    |                |
| SPW Ground stud        | RF common ground stud                        | 85.6 (2200)    |                |
| PGR, top panel         | GOC, rear panel                              | 1133.9 (28800) |                |
| PGR, top panel         | HEC, top panel                               | 346.5 (8800)   |                |
| PGR, top panel         | Magnet room door switch                      | 960.6 (24400)  |                |
| PGR, top panel         | PEN cabinet, top panel                       | 370.1 (9400)   | 645.7 (16400)  |
| PGR, top panel         | SPW, bottom edge                             | 385.8 (9800)   | 661.4 (16800)  |
| PGR, top panel         | RF common ground stud                        | 362.2 (9200)   | 637.8 (16200)  |
| SPW, bottom edge       | CRY                                          | 551.2 (14000)  |                |
| SPW, bottom edge       | E-off switch, Control room or Equipment room | 385.8 (9800)   | 661.4 (16800)  |
| SPW, bottom edge       | GOC, rear panel                              | 1149.6 (29200) |                |
| HEC, top panel         | SPW, bottom edge                             | 244 (6200)     | 677 (17200)    |
| <b>Magnet Room</b>     |                                              |                |                |
| PEN panel, bottom edge | MAG, gradient cable clamp block              | 354.3 (9000)   | 511.8 (13000)  |
| SPW, bottom edge       | MAG, gradient cable clamp block              | 354.3 (9000)   | 511.8 (13000)  |
| SPW, bottom edge       | E-off switch, Magnet Room                    | 1173.2 (29800) |                |
| GOC, rear panel        | MAG, gradient cable clamp block              | 1496.1 (38000) | 1653.5 (42000) |
| MRU                    | MAG, gradient cable clamp block              | 972.4 (24700)  |                |
| RF common ground stud  | MAG, gradient cable clamp block              | 393.7 (10000)  | 551.2 (14000)  |

### 1.3 Brainwave Option

Table 6-4: Brainwave Option Usable Cable Lengths

| Cable Part number | Point A            | Point B            | Site Option A:<br>Short ER, Short SR | Site Option B: Long<br>ER, Short SR | Site Option C:<br>Short ER, Long SR |
|-------------------|--------------------|--------------------|--------------------------------------|-------------------------------------|-------------------------------------|
|                   |                    |                    | in. (mm)                             |                                     |                                     |
|                   | PEN Wall/Cabinet   | Brainwavew Cabinet |                                      |                                     |                                     |
| LG408             | Empty 9 pin D-sub  | Lumina Controller  | 720 (18,288)                         |                                     |                                     |
| LU001-60          | J16 in PEN Cabinet | Lumina Controller  | 720 (18,288)                         |                                     |                                     |

### 1.4 Magnetic Resonance Elastography (MRE) Option

Table 6-5: MRE Option Usable Cable Lengths

| Cable Part number | Point A | Point B | Site Option A:<br>Short ER, Short SR | Site Option B: Long<br>ER, Short SR | Site Option C:<br>Short ER, Long SR |
|-------------------|---------|---------|--------------------------------------|-------------------------------------|-------------------------------------|
|                   |         |         | in. (mm)                             |                                     |                                     |

|                       |                          |                          |                                           |
|-----------------------|--------------------------|--------------------------|-------------------------------------------|
| 1 in (2.54 mm) Tubing | Resoundant Active Driver | Magnet (Isoceter)        | Nominal 288 (7315)<br>Maximum 396 (10058) |
| BNC                   | Resoundant Active Driver | PEN Cabinet J6           | 600 (15,240)                              |
| Ethernet              | Resoundant Active Driver | Ethernet Hub in PGR      | 600 (15,240)                              |
| Power                 | Resoundant Active Driver | Customer Supplied Outlet | 60 Hz: 240 (6096)<br>50 Hz: 300 (7620)    |

## 2 MR System Interconnects Routing Requirements


### 2.1 General Requirements

1. The customer is responsible for the purchase and installation of any Interconnect Cable Routing Mechanism (ICRM)
2. Any ICRM may be used provided it meets all MR system requirements and any applicable local and national codes
3. The ICRM must be nonferrous (e.g., composites or aluminum)
4. The ICRM must meet the minimum cable bend radius. See [MR System Cable Specifications](#)
5. Each ICRM must support a weight of at least 50 lbs/ft (74.8 kg/m)

**NOTE:** This weight requirement does not include any force applied during cable installation in addition to the cable weight.

6. Two ICRMs must be used, each at least 18 in. (457 mm) wide
7. Power and signal cables must be routed in an ICRM separate from the air, water, and, cryogen line ICRM. See [Illustration 6-1](#)
8. The gradient and the RF common ground cables must be routed in a separate subsection of the ICRM
9. The following cable groups must be routed in separate sections of the ICRM:
  - a. Coax, RF, and AC power cables with jacket rating of 600V and above
  - b. Data and fiber optic and 300V coax or RF clock cables
  - c. 300V signal, 300V power, and 300V power/signal cables

**Illustration 6-1: Cable Groupings (Magnet and Equipment Room)**

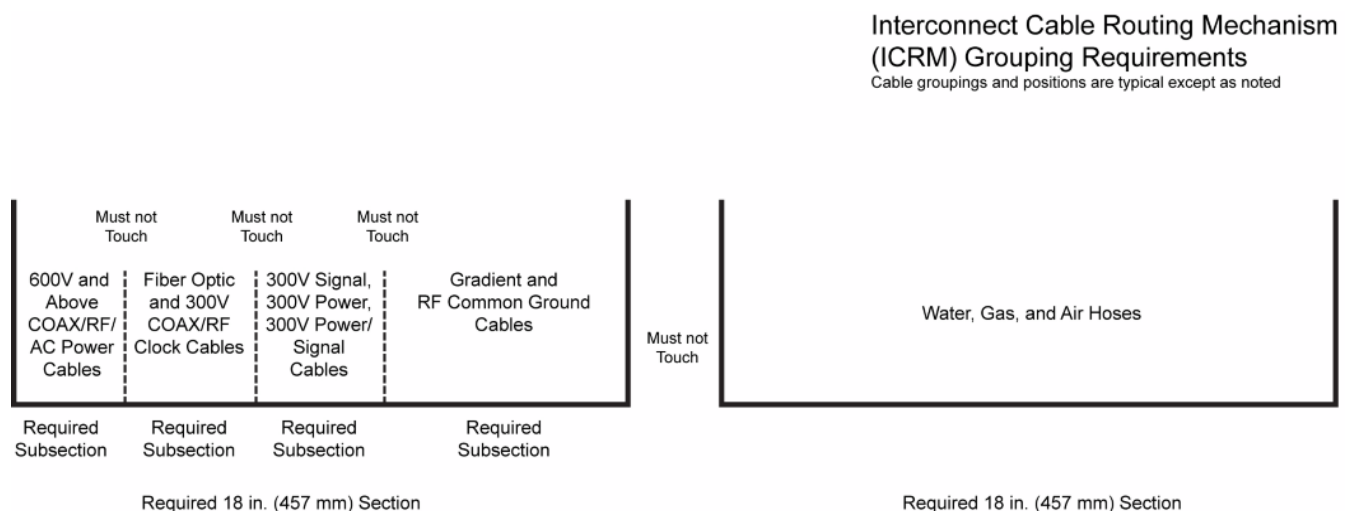

10. For multiple MR system installations, cables from different MR systems must not share ICRMs

### 2.2 Magnet Room Requirements

1. ICRMs must attach to the RF Shield Room and installation must be coordinated with the RF shield vendor


2. Cables must not be routed over the magnet (see [Chapter 2, MR Suite Minimum Room Size Requirements](#) for ceiling exclusion area)
3. The magnet end of the ICRM must end 44.45 in. ( $\pm 0.5$ ) [1205 mm ( $\pm 12$ )] from geometric isocenter. See [Illustration 6-2](#)
4. ICRMs must not be stacked (side-to-side only)
5. ICRMs must not touch each other
6. Ceiling grid work, medical gas lines, lighting fixtures, etc. must not touch MR system cabling or ICRMs and must not interfere with the service area around the ICRM.
7. Cables must be accessible for maintenance at all points along the route

**NOTICE**

The minimum ceiling height dimensions listed below are for ICRM routing only. Overall ceiling height must also take into account minimum service area requirements [MR Suite Minimum Room Size Requirements](#) and cryogen venting requirements [Chapter 3, Magnet Room Venting Requirements](#).

- a. The top of the ICRM must be between 106 in. (2692 mm) and 116 in. (2946 mm) above the finished floor
- b. One-side accessible only ICRMs require a minimum of 14 in. (355 mm) from the top of the mechanism to any object above the ICRM and the ICRMs must be at least 6 in. (152 mm) apart
- c. ICRMs must be accessible on both sides require a minimum of 10 in. (254 mm) from the top of the ICRM to any object above the ICRM


Illustration 6-2: Interconnect Cable Routing Mechanism Positioning

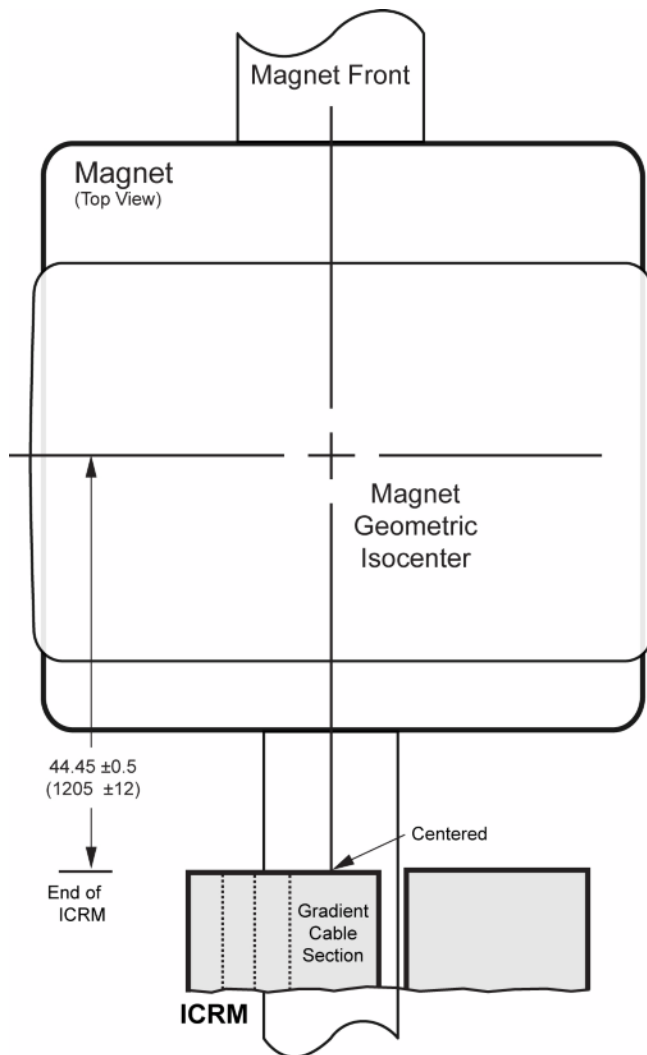

ICRM Magnet-end Positioning

## Interconnect Cable Routing Mechanism (ICRM) Requirements

All Dimensions are in inches  
Bracketed dimensions are in millimeters

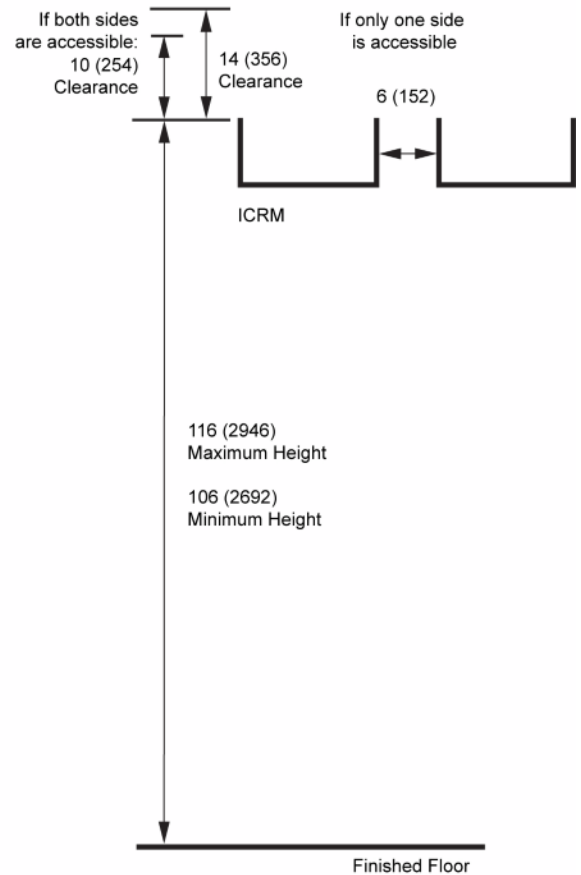

ICRM Height

8. All Magnet Room interconnects must route overhead and enter the Magnet enclosure from the top
9. All cables must approach the magnet from the rear wall of the scan room
10. The magnet-end subsection of the ICRM routing the gradient cables must be aligned to the center of the magnet (when viewed from the top)
11. All electrical and mechanical connections and fasteners must be tightened and secured to supplier specifications to prevent broadband interference
12. Excess cable length in the Magnet room must be stored either in the:
  - a. Penetration Panel closet
  - b. Magnet Room ICRM (excess cable must be at least 36 in. (914 mm) from the magnet end of the ICRM)



## 2.3 Equipment Room Requirements

1. All equipment interconnects must route overhead
2. ICRMs must route directly above Equipment room components and cables/hoses must drop directly to the top of the cabinets (See individual components in Chapter 4 for ICRM height requirements)
3. ICRMs do not have height restrictions except above Equipment room components (however, all service access requirements must be met)
4. If the ICRMs are routed above a finished ceiling, the customer must provide access for the cable drops
5. If ICRMs are routed above a finished ceiling in the Equipment room, there must be at least 24 in. (610 mm) of clear space above the finished ceiling
6. Excess cable length in the equipment must be stored in the Equipment room

### 3 MR System Cable Specifications

MR450w cables are marked with an X in the column titled 450w. The other cables are not used with the MR450w system.

Illustration 6-3: System Cable List (Page 1)

|     |     |      |             | Connector                                                                                    |       |    | Cable |    |     | Bend |    |    | Cond-<br>uctors | AWG | Certi-<br>fication | Volt<br>Rating<br>(volts) | Temp<br>Rating<br>(C) | Notes                            |
|-----|-----|------|-------------|----------------------------------------------------------------------------------------------|-------|----|-------|----|-----|------|----|----|-----------------|-----|--------------------|---------------------------|-----------------------|----------------------------------|
| 450 | 750 | 450w | Part Number | Description                                                                                  | in    | mm | in    | mm | in  | mm   | mm | mm |                 |     |                    |                           |                       |                                  |
| X   | X   | X    | 2200835     | Run E3009, HEC to CRY, INPUT POWER, SUMITOMO CSW-71D, 10000mm                                | 1.5   | 38 | 0.82  | 21 | 7.0 | 178  |    |    | 4               | 12  | LFNC               | 600                       | 105°                  | Four 12AWG wires in conduit      |
| X   | X   | X    | 2155316     | Run 623, CRY to SPW, COLD HEAD CONTROL, LOC, 15000mm                                         | 1.25  | 32 | 0.28  | 7  | 7.0 | 178  |    |    | 4               | 16  | CM, FT4            | 300                       | 90°                   |                                  |
| X   | X   | X    | 2218292     | Run 624, SPW to MAG, COLD HEAD CONTROL, LOC, 15000mm                                         | 1.25  | 32 | 0.28  | 7  | 7.0 | 178  |    |    | 4               | 16  | CM, FT4            | 300                       | 90°                   |                                  |
| X   | X   | X    | 5191457     | Run M3030, SPW-J115 to E-JOFF, 30000mm                                                       | 1.5   | 38 | 0.305 | 8  | 6.0 | 152  |    |    | 9               | 22  | CL3, FT4           | 300                       | 80°                   |                                  |
| X   | X   | X    | 2384104-2   | Run 1152, PPI-45-P1, TO MAG2A12A3, 15240mm                                                   | 1.75  | 44 | 0.58  | 15 | 9.0 | 229  |    |    | 13              | 16  | CL3, FT4           | 300                       | 105°                  | Higher order shim cable (option) |
| X   | X   | X    | 5168522-2   | Run E3008, PGR-J4 to PEN-J9, CAN, 10000mm                                                    | 1.5   | 38 | 0.305 | 8  | 6.0 | 152  |    |    | 9               | 22  | CL3, FT4           | 300                       | 80°                   |                                  |
| X   | X   | X    | 5168522-3   | Run E3008, PGR-J9 to PEN-J9, CAN, 17000mm                                                    | 1.5   | 38 | 0.305 | 8  | 6.0 | 152  |    |    | 9               | 22  | CL3, FT4           | 300                       | 80°                   |                                  |
| X   | X   | X    | 5168749-2   | Run M1302, PEN-J29 to MAG-TDM-J11, Receiver Clock, Twimax, 10000mm                           | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5168749-3   | Run M1304, PEN-J30 to MAG-J8, Rx LO - Exciter Reference Clock, Twimax, 10000mm               | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5168749-6   | Run M1302, PEN-J29 to MAG-TDM-J11, Receiver Clock, Twimax, 14000mm                           | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5168749-7   | Run M1304, PEN-J30 to MAG-J8, Rx LO - Exciter Reference Clock, Twimax, 14000mm               | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5330926-2   | Run M1302, PEN-J29 to MAG-TDM-J11, Receiver Clock, Twimax, 11000mm                           | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5330926-6   | Run M1302, PEN-J29 to MAG-TDM-J11, Receiver Clock, Twimax, 15000mm                           | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5330926-3   | Run M1304, PEN-J30 to PED-HUB-SL0T11-J17, Rx LO - Exciter Reference Clock, Twimax, 13500mm   | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5330926-7   | Run M1304, PEN-J30 to PED-HUB-SL0T11-J17, Rx LO - Exciter Reference Clock, Twimax, 17500mm   | 1     | 25 | 0.33  | 8  | 5.0 | 127  |    |    | 2               | 20  | CL2, FT4           | 300                       | 75°                   | twimaxial cable                  |
| X   | X   | X    | 5167907-2   | Run E3011, PGR-J2 to OM1-Host, Ethernet, 30000mm                                             | 0.625 | 16 | 0.25  | 6  | 5.0 | 127  |    |    | 8               | 26  | CMP                | 300                       | 75°                   |                                  |
| X   | X   | X    | 5167907-4   | Run E3013, PGR-J6 to HEC-CAT 5, Cooling Control Ethernet, 10000mm                            | 0.625 | 16 | 0.25  | 6  | 5.0 | 127  |    |    | 8               | 26  | CMP                | 300                       | 75°                   |                                  |
| X   | X   | X    | 5167907-5   | Run E3034 MON-J5 to customer-supplied Network, 30000mm                                       | 0.625 | 16 | 0.25  | 6  | 5.0 | 127  |    |    | 8               | 26  | CMP                | 300                       | 75°                   |                                  |
| X   | X   | X    | 5168163     | Run E1001, PGR-J19 to SPW-J106, 1.5T RF Transmit, 10327mm ref                                | 2.85  | 72 | 0.9   | 23 | 9.0 | 229  |    |    | 1               | N/A | CMR, FT4           | 5000                      | 85°                   | LMR-900 coaxial cable            |
| X   | X   | X    | 5168163-2   | Run M1001, SPW-J106 to RUN M1319, 1.5T RF Transmit, 10327 mm ref                             | 2.85  | 72 | 0.9   | 23 | 9.0 | 229  |    |    | 1               | N/A | CMR, FT4           | 5000                      | 85°                   | LMR-900 coaxial cable            |
| X   | X   | X    | 5168163-3   | Run M1001, SPW-J106 to RUN M1319, 1.5T RF Transmit, 14457mm ref                              | 2.85  | 72 | 0.9   | 23 | 9.0 | 229  |    |    | 1               | N/A | CMR, FT4           | 5000                      | 85°                   | LMR-900 coaxial cable            |
| X   | X   | X    | 5168163-4   | Run E1001, PGR-J19 to SPW-J106, 1.5T RF Transmit, 16523mm ref                                | 2.85  | 72 | 0.9   | 23 | 9.0 | 229  |    |    | 1               | N/A | CMR, FT4           | 5000                      | 85°                   | LMR-900 coaxial cable            |
| X   | X   | X    | 5344701     | Run M1001, SPW-J106 to MAG-HYB-J1, 1.5T RF Body Transmit, 12394mm ref, Aluminum braid shield | 2.85  | 72 | 0.9   | 23 | 9.0 | 229  |    |    | 1               | N/A | CMR, FT4           | 5000                      | 85°                   | LMR-900 coaxial cable            |
| X   | X   | X    | 5344701-2   | Run M1001, SPW-J106 to MAG-HYB-J1, 1.5T RF Body Transmit, 16525mm ref, Aluminum braid shield | 2.85  | 72 | 0.9   | 23 | 9.0 | 229  |    |    | 1               | N/A | CMR, FT4           | 5000                      | 85°                   | LMR-900 coaxial cable            |
| X   | X   | X    | 5330729-5   | Run E1001, PGR-J19 to SPW-J106, 1.5T RF Transmit, 10327mm ref, Aluminum braid shield         | 2.85  | 72 | 0.9   | 23 | 9.0 | 229  |    |    | 1               | N/A | CMR, FT4           | 5000                      | 85°                   | LMR-900 coaxial cable            |

Illustration 6-4: System Cable List (Page 2)

| 450 | 750 | 450w | Part Number | Description                                                                          | Connector | Cable |       |    | Bend |     | Cond-uctors | AWG                   | Certi-fication | Volt Rating (volts) | Temp Rating (C) | Notes                 |
|-----|-----|------|-------------|--------------------------------------------------------------------------------------|-----------|-------|-------|----|------|-----|-------------|-----------------------|----------------|---------------------|-----------------|-----------------------|
|     |     |      |             |                                                                                      | In        | mm    | In    | mm | In   | mm  |             |                       |                |                     |                 |                       |
|     |     |      | 5330729-6   | Run E1001, PGR-J19 to SPW-J106, 1.5 RF Transmitt, 1692mm ref. Aluminum braid shield. | 2.85      | 72    | 0.9   | 23 | 9.0  | 229 | 1           | N/A                   | CMR, FT4       | 5000                | 85°             | LMR-800 coaxial cable |
| X   | X   | X    | 5168733     | Run E0003, PGR-J102 to OW-PDM-J1, Host Power, 30000mm                                | 1.4       | 36    | 0.45  | 11 | 6.0  | 152 | 3           | 10                    | UL AWM         | 600                 | 90°             |                       |
| X   | X   | X    | 5168845     | Run E3028, HEC-PLC-RS232 to PGR-J3, Cooling Control Communication, 10000mm           | 0.625     | 16    | 0.25  | 6  | 6.0  | 152 | 8           | 26                    | CM, FT4        | 300                 | 75°             |                       |
| X   | X   | X    | 5168970     | Run E3030, SPW-J115 to E-Off SW, E-Off 10000mm                                       | 1.5       | 38    | 0.3   | 8  | 5.0  | 127 | 2           | 22                    | CMG, FT4       | 300                 | 105°            |                       |
| X   | X   |      | 5168970-2   | Run E3003, PGR-J7 to MDP-TB1-14-17, E-Off 10000mm                                    | 1.5       | 38    | 0.3   | 8  | 5.0  | 127 | 2           | 22                    | CMG, FT4       | 300                 | 105°            |                       |
| X   | X   |      | 5168970-3   | Run E3027, HEC-MDP to MDP-TB1-17-18, E-Off 10000mm                                   | 1.5       | 38    | 0.3   | 8  | 5.0  | 127 | 2           | 22                    | CMG, FT4       | 300                 | 105°            |                       |
| X   | X   | X    | 5168970-4   | Run E3014, PGR-J20 to RF Door SW, E-Stop, 25000mm                                    | 1.5       | 38    | 0.3   | 8  | 5.0  | 127 | 2           | 22                    | CMG, FT4       | 300                 | 105°            |                       |
| X   | X   | X    | 5168970-5   | Run E3030, E-Off Switch to SPW-J115, 17000mm                                         | 1.5       | 38    | 0.3   | 8  | 5.0  | 127 | 2           | 22                    | CMG, FT4       | 300                 | 105°            |                       |
| X   | X   | X    | 5169207-2   | Run E0004, PGR-J16 to PEN-J1, PEN Power, 10000mm                                     | 1.75      | 44    | 0.57  | 14 | 6.0  | 152 | 5           | 10                    | UL AWM         | 600                 | 90°             |                       |
| X   | X   | X    | 5169207-4   | Run E0004, PGR-J16 to Pen-J1, Pen Power, 17000mm                                     | 1.75      | 44    | 0.57  | 14 | 6.0  | 152 | 5           | 10                    | UL AWM         | 600                 | 90°             |                       |
| X   | X   | X    | 5169227-10  | Run M3031, SPW-J118 to MAG-J1, Console Intercom, 10000mm                             | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22                    | CL3, FT4       | 300                 | 80°             |                       |
| X   | X   | X    | 5169227-11  | Run M3031, SPW-J118 to MAG-J1, Console Intercom, 14000mm                             | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22                    | CL3, FT4       | 300                 | 80°             |                       |
| X   | X   | X    | 5169227-2   | Run E3031, SPW-J118 to OW-J7, Patient Intercom, 30000mm                              | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22                    | CL3, FT4       | 300                 | 80°             |                       |
| X   | X   | X    | 5169227-5   | Run M3312, PEN-J41 to MAG-J6, PAC Power, 10000mm                                     | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22                    | CL3, FT4       | 300                 | 80°             |                       |
| X   | X   | X    | 5169227-6   | Run M3311, Mag-J16 to Pen-J36, Long Drive Power, 10000mm                             | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22                    | CL3, FT4       | 300                 | 80°             |                       |
| X   | X   | X    | 5169227-7   | Run M3312, PEN-J41 to MAG-J6, PAC Power, 14000mm                                     | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22                    | CL3, FT4       | 300                 | 80°             |                       |
| X   | X   |      | 5169227-8   | Run M3311, Pen-J36 to Mag-J16, Long Drive Power, 14000mm                             | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22                    | CL3, FT4       | 300                 | 80°             |                       |
|     |     | X    | 5329980     | Run M3311, PEN-J36 to PED-MTR-J1, Long Drive Power, 15500mm                          | 2.5       | 64    | 0.39  | 10 | 3.9  | 99  | 6           | 4/C 14AWG + 2/C 22AWG | CMG, FT4       | 300                 | 80°             |                       |
| X   | X   | X    | 5329990-2   | Run M3311, PEN-J36 to PED-MTR-J1, Long Drive Power, 19500mm                          | 2.5       | 64    | 0.39  | 10 | 3.9  | 99  | 6           | 4/C 14AWG + 2/C 22AWG | CMG, FT4       | 300                 | 80°             |                       |
| X   | X   | X    | 5169277     | Run E3006, PGR-J5 to PEN-J14, HART, 10000mm                                          | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   | X    | 5169277-10  | Run M3316, Pen-J54 to Mag-J7, 14000mm                                                | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   | X    | 5169277-11  | Run M3315, Pen-J54 to Mag-J17, 14000mm                                               | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   | X    | 5169277-12  | Run M3022, SPW-J116 to MAG-P403, 14000 mm                                            | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   | X    | 5169277-2   | Run E3015, PGR-J17 to OW-F-J6, E-Stop, 30000mm                                       | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   | X    | 5169277-4   | Run E3017, PEN-J10 to OW-F-J5, E-Stop, 30000mm                                       | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   |      | 5169277-6   | Run M3316, PEN-J34 to MAG-J7, TRDD Unblank, 10000mm                                  | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   |      | 5169277-7   | Run M3315, PEN-J54 to MAG-J17, Half's Signal, 10000mm                                | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |
| X   | X   | X    | 5169277-9   | Run E3006, PGR-J5 to Pen-J14, 17000 mm                                               | 1.5       | 38    | 0.305 | 8  | 6.0  | 152 | 9           | 22                    | CL3, FT4       | 300                 | 90°             |                       |

Illustration 6-5: System Cable List (Page 3)

| 450 | 750 | 450w | Part Number | Description                                                        | Connector | Cable | Bend  | Condi-<br>uctors | AWG | Certi-<br>fication | Volt<br>Rating<br>(volts) | Temp<br>Rating<br>(C) | Notes                 |
|-----|-----|------|-------------|--------------------------------------------------------------------|-----------|-------|-------|------------------|-----|--------------------|---------------------------|-----------------------|-----------------------|
|     |     |      |             |                                                                    | In        | mm    | In    | mm               | In  | mm                 |                           |                       |                       |
|     |     |      |             |                                                                    | 1.5       | 38    | 0.305 | 8                | 6.0 | 152                | 9                         | 90°                   |                       |
|     |     |      | 5330927-6   | Run M3316, PEN-J34 to PED-HUB-SLOT11-J10, TR DD Unibank, 13500mm.  |           |       |       |                  |     |                    |                           |                       |                       |
|     |     |      | 5330927-10  | Run M3316, PEN-J34 to PED-HUB-SLOT11-J10, 17500mm.                 | 1.5       | 38    | 0.305 | 8                | 6.0 | 152                | 9                         | 90°                   |                       |
|     |     |      | 5330927-7   | Run M3315, PEN-J54 to PED-TIF-J4, Hall's Signal, 14500mm.          | 1.5       | 38    | 0.305 | 8                | 6.0 | 152                | 9                         | 90°                   |                       |
|     |     |      | 5330927-11  | Run M3315, PEN-J54 to PED-TIF-J4, 18500mm.                         | 1.5       | 38    | 0.305 | 8                | 6.0 | 152                | 9                         | 90°                   |                       |
|     |     |      | 5169388     | Run E3009 PGR-J48 to Pen-J13, 10000 mm                             | 3         | 76    | 0.52  | 13               | 6.0 | 152                | 37                        | 80°                   |                       |
|     |     |      | 5169388-2   | Run E3009 PGR-J48 to Pen-J13, 17000 mm                             | 3         | 76    | 0.52  | 13               | 6.0 | 152                | 37                        | 80°                   |                       |
|     |     |      | 5169404     | Run E3018, PEN-J12 to OW-F17, HART-RS422, 30000mm                  | 1.75      | 44    | 0.34  | 9                | 6.0 | 152                | 15                        | 80°                   |                       |
|     |     |      | 5169404-2   | Run E3002, PGR-J103 to PEN-J11, E-Stop, 10000mm                    | 1.75      | 44    | 0.34  | 9                | 6.0 | 152                | 15                        | 80°                   |                       |
|     |     |      | 5169404-3   | Run M3314, PEN-J42 to MAG-J2, LCD Power, 10000mm                   | 1.75      | 44    | 0.34  | 9                | 6.0 | 152                | 15                        | 80°                   |                       |
|     |     |      | 5169404-4   | Run M3314, Pen-J42 to MAG-J2, 14000 mm                             | 1.75      | 44    | 0.34  | 9                | 6.0 | 152                | 15                        | 80°                   |                       |
|     |     |      | 5169404-5   | Run E3002 PGR-J103 to Pen-J11, 17000 mm                            | 1.75      | 44    | 0.34  | 9                | 6.0 | 152                | 15                        | 80°                   |                       |
|     |     |      | 5169404-7   | Run E3020 Pen-J13 to Men-J10, 17000mm                              | 1.75      | 44    | 0.34  | 9                | 6.0 | 152                | 15                        | 80°                   |                       |
|     |     |      | 5169588     | Run M3310 PEN-J45 to PED-HUB-SLOT12-J4, RF Power, 13500mm          | 3.5       | 89    | 0.44  | 11               | 6.8 | 173                | 25                        | 80°                   |                       |
|     |     |      | 5169588-2   | Run M3310 PEN-J45 to PED-HUB-SLOT12-J4, RF Power, 17500mm          | 3.5       | 89    | 0.44  | 11               | 6.8 | 173                | 25                        | 80°                   |                       |
|     |     |      | 5169804-11  | Run E1006, PGR-J9 to SPW-J107, Dummy Load SAR Monitoring, 10000mm  | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-12  | Run E1015, PGR-J10 to SPW-J108, Dummy Load SAR Monitoring, 10000mm | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-13  | Run M1303, PEN-J32 to MAG-J13, NB Exciter, 10000mm                 | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-20  | Run M1006, SPW-J107 to PED-DUM-J13, Dummyload, 3T, 13500mm         | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-21  | Run M1015, SPW-J108 to PED-DUM-J14, Dummyload, 3T, 13500mm         | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-25  | Run E1005, PGR-J11 to PEN-J15, Tx Output, 17000mm                  | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-26  | Run E1006, PGR-J9 to SPW-J107, Dummy Load SAR Monitoring, 17000mm  | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-27  | Run E1015, PGR-J10 to SPW-J108, Dummy Load SAR Monitoring, 17000mm | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-28  | Run M1303, PEN-J32 to MAG-J13, NB Exciter, 14000mm                 | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-29  | Run M1006, SPW-J107 to PED-DUM-J13, Dummyload, 3T, 17500mm         | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-3   | Run E1005, PGR-J11 to PEN-J15, Tx Output, 10000mm                  | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5169804-30  | Run M1015, SPW-J108 to PED-DUM-J14, Dummyload, 3T, 17500mm         | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5329514-13  | Run M1303, PEN-J32 to PED-HUB-SLOT11-J16, NB Exciter, 14500mm      | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |
|     |     |      | 5329514-28  | Run M1303, PEN-J32 to PED-HUB-SLOT11-J16, NB Exciter, 18500mm      | 0.625     | 16    | 0.195 | 5                | 2.0 | 51                 | 1                         | 100°                  | LMR-195 coaxial cable |

### Illustration 6-6: System Cable List (Page 4)

[illegible]



### Illustration 6-7: System Cable List (Page 5)

[illegible]

### Illustration 6-8: System Cable List (Page 6)

|     |     |      |             |                                                                                               | Connector                                              | Cable | Braid | Cond-<br>uctors | AWG | Certifi-<br>cation | Volt<br>Rating<br>(volts) | Temp<br>Rating<br>(C) | Notes                  |                    |     |                       |                                                        |
|-----|-----|------|-------------|-----------------------------------------------------------------------------------------------|--------------------------------------------------------|-------|-------|-----------------|-----|--------------------|---------------------------|-----------------------|------------------------|--------------------|-----|-----------------------|--------------------------------------------------------|
| 450 | 750 | 450w | Part Number | Description                                                                                   | In                                                     | mm    | In    | mm              | In  | mm                 |                           |                       |                        |                    |     |                       |                                                        |
|     |     |      |             | Run M1305, SPW-J122 to PED-H-1, P1 port RF Transmit, 14306mm reference, Aluminum braid shield | 1                                                      | 25    | 0.59  | 15              | 6.0 | 152                | 1                         | N/A                   | CMR, FT4               | 4000               | 85° | LMR-600 coaxial cable |                                                        |
|     |     |      | 5329340-8   | Run M1305, SPW-J122 to PED-H-1, P1 port RF Transmit, 18393mm reference, Aluminum braid shield | 1                                                      | 25    | 0.59  | 15              | 6.0 | 152                | 1                         | N/A                   | CMR, FT4               | 4000               | 85° | LMR-600 coaxial cable |                                                        |
|     |     |      | 5329340-9   | Run M1306, SPW-J120 to PED-H-2, A port RF Transmit, 13865mm reference, Aluminum braid shield  | 1                                                      | 25    | 0.59  | 15              | 6.0 | 152                | 1                         | N/A                   | CMR, FT4               | 4000               | 85° | LMR-600 coaxial cable |                                                        |
|     |     |      | 5329340-10  | Run M1306, SPW-J120 to PED-H-2, A port RF Transmit, 18052mm reference, Aluminum braid shield  | 1                                                      | 25    | 0.59  | 15              | 6.0 | 152                | 1                         | N/A                   | CMR, FT4               | 4000               | 85° | LMR-600 coaxial cable |                                                        |
|     |     |      | 5329340-11  | Run E1002 SPW-J128 to PGR-J18, RF Transmit, 10216mm, Aluminum braid shield                    | 1                                                      | 25    | 0.59  | 15              | 6.0 | 152                | 1                         | N/A                   | CMR, FT4               | 4000               | 85° | LMR-600 coaxial cable |                                                        |
|     |     |      | 5329340-12  | Run E1002 SPW-J128 to PGR-J18, RF Transmit, 16348mm, Aluminum braid shield                    | 1                                                      | 25    | 0.59  | 15              | 6.0 | 152                | 1                         | N/A                   | CMR, FT4               | 4000               | 85° | LMR-600 coaxial cable |                                                        |
|     |     |      | X           | Run EX305, CRY to MON-FJ4, Flow Monitor-Sunittomo, 10000mm                                    | 1.75                                                   | 44    | 0.35  | 9               | 6.0 | 152                | 10                        | 22                    | CMG, FT4               | 300                | 80° | LMR-600 coaxial cable |                                                        |
|     |     |      | X           | 5174110                                                                                       |                                                        |       |       |                 |     |                    |                           |                       |                        |                    |     |                       |                                                        |
|     |     |      | X           | 5174662                                                                                       | Run M2025 Mag to Pen, 12500mm                          | 2.4   | 61    | 0.82            | 21  | 7.0                | 178                       | 8                     | N/A                    | LFNC               | N/A | 85°                   | optical fibers in conduit, 4 duplex connectors bundled |
|     |     |      | X           | 5174662-2                                                                                     | Run M2025 Mag to Pen, 16500mm                          | 2.4   | 61    | 0.82            | 21  | 7.0                | 178                       | 8                     | N/A                    | LFNC               | N/A | 85°                   | optical fibers in conduit, 4 duplex connectors bundled |
|     |     |      | X           | 5344936                                                                                       | Run M2025 MAG to PEN, CAN fiber optic, 12500mm         | 2.4   | 61    | 0.82            | 21  | 7.0                | 178                       | 8                     | N/A                    | LFNC               | N/A | 85°                   | optical fibers in conduit, 4 duplex connectors bundled |
|     |     |      | X           | 5344936-2                                                                                     | Run M2025 MAG to PEN, CAN fiber optic, 16500mm         | 2.4   | 61    | 0.82            | 21  | 7.0                | 178                       | 8                     | N/A                    | LFNC               | N/A | 85°                   | optical fibers in conduit, 4 duplex connectors bundled |
|     |     |      | X           | 5194585                                                                                       | Run M3317 Gradient Cable X, Scan Room, Short, 10000mm  | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194585-2                                                                                     | Run M3318 Gradient Cable Y, Scan Room, Short, 10000mm  | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194585-3                                                                                     | Run M3319 Gradient Cable Z, Scan Room, Short, 10000mm  | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194585-4                                                                                     | Run M3317 Gradient Cable X, Scan Room, Long, 14000mm   | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194585-5                                                                                     | Run M3318 Gradient Cable Y, Scan Room, Long, 14000mm   | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194585-6                                                                                     | Run M3319 Gradient Cable Z, Scan Room, Long, 14000mm   | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194889                                                                                       | Run E3317 Gradient Cable X, Equip Room, Short, 10000mm | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194889-2                                                                                     | Run E3318 Gradient Cable Y, Equip Room, Short, 10000mm | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194889-3                                                                                     | Run E3317 Gradient Cable X, Equip Room, Short, 10000mm | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194889-4                                                                                     | Run E3318 Gradient Cable Y, Equip Room, Long, 17000mm  | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194889-5                                                                                     | Run E3319 Gradient Cable Z, Equip Room, Long, 17000mm  | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |
|     |     |      | X           | 5194889-6                                                                                     | Run E3319 Gradient Cable Z, Equip Room, Long, 17000mm  | 2     | 51    | 2               | 51  | 13.0               | 330                       | 5                     | 4/C 20 AWG + 1/C 6 AWG | UL Type TC-ER, FT4 | 600 | 90°                   |                                                        |



Illustration 6-9: System Cable List (Page 7)

| 450 | 750 | 450w | Part Number | Description                                         | Connector |    | Cable |    | Bend |     | Conductors | AWG | Certification    | Volt Rating (volts) | Temp Rating (C) | Notes                                                       |
|-----|-----|------|-------------|-----------------------------------------------------|-----------|----|-------|----|------|-----|------------|-----|------------------|---------------------|-----------------|-------------------------------------------------------------|
|     |     |      |             |                                                     | In        | mm | In    | mm | In   | mm  |            |     |                  |                     |                 |                                                             |
| X   | X   |      | 5189365     | Run E1002 SPW-J128 to PGR-J18, RF Transmit, 10218mm | 1.36      | 35 | 0.9   | 23 | 9.0  | 229 | 1          | N/A | CMR, FT4         | 5000                | 85°             | LMR-900 coaxial cable                                       |
| X   | X   |      | 5189365-2   | Run E1002 SPW-J128 to PGR-J18, RF Transmit, 16348mm | 1.36      | 35 | 0.9   | 23 | 9.0  | 229 | 1          | N/A | CMR, FT4         | 5000                | 85°             | LMR-900 coaxial cable                                       |
| X   | X   | X    | 5191204     | Run E2002 PGR-J1 to PEN-J19, 10000mm                | 1.2       | 30 | 0.32  | 8  | 4.8  | 122 | 2          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 2 simplex connectors bundled on each end    |
| X   | X   | X    | 5191204-2   | Run E2002 PGR-J1 to PEN-J19, 17000mm                | 1.2       | 30 | 0.32  | 8  | 4.8  | 122 | 2          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 2 simplex connectors bundled on each end    |
| X   | X   | X    | 5191204-3   | Run E2023 PGR-J21 to PEN-J18, 10000mm               | 1.2       | 30 | 0.32  | 8  | 4.8  | 122 | 2          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 2 simplex connectors bundled on each end    |
| X   | X   | X    | 5191204-4   | Run E2023 PGR-J21 to PEN-J18, 17000mm               | 1.2       | 30 | 0.32  | 8  | 4.8  | 122 | 2          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 2 simplex connectors bundled on each end    |
| X   | X   | X    | 5191301-2   | Run P2003 PGR-J15 to Mag-J20, 31000mm               | 1.2       | 30 | 0.38  | 10 | 5.7  | 145 | 6          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 3 duplex connectors bundled on each end     |
| X   | X   | X    | 5191301-3   | Run P2004 PGR-J14 to Mag-J21, 20000mm               | 1.2       | 30 | 0.38  | 10 | 5.7  | 145 | 6          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 3 duplex connectors bundled on each end     |
| X   | X   | X    | 5191301-4   | Run P2004 PGR-J14 to Mag-J21, 31000mm               | 1.2       | 30 | 0.38  | 10 | 5.7  | 145 | 6          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 3 duplex connectors bundled on each end     |
| X   | X   | X    | 5191924     | Run E1003, PGR-J12 to PEN-J5, Body TR, 10000mm      | 0.625     | 16 | 0.24  | 6  | 2.5  | 64  | 1          | N/A | CMR, FT4         | 1500                | 85°             | LMR-240 coaxial cable                                       |
| X   | X   | X    | 5191924-2   | Run E1004 PGR-J13 to Pen-J6 Head TR, 10000mm        | 0.625     | 16 | 0.24  | 6  | 2.5  | 64  | 1          | N/A | CMR, FT4         | 1500                | 85°             | LMR-240 coaxial cable                                       |
| X   | X   | X    | 5191924-3   | Run E1003 PGR-J12 to PEN-J5 Body TR, 17000mm        | 0.625     | 16 | 0.24  | 6  | 2.5  | 64  | 1          | N/A | CMR, FT4         | 1500                | 85°             | LMR-240 coaxial cable                                       |
| X   | X   | X    | 5191924-4   | Run E1004, PGR-J13 to PEN-J6, Head TR, 17000mm      | 0.625     | 16 | 0.24  | 6  | 2.5  | 64  | 1          | N/A | CMR, FT4         | 1500                | 85°             | LMR-240 coaxial cable                                       |
| X   | X   | X    | 5192963     | Run E2020 PGR-J22 to PEN-J17, 10000mm               | 2.4       | 61 | 0.82  | 21 | 7.0  | 178 | 6          | N/A | LFNC             | N/A                 | 85°             | optical fibers in conduit, 3 connectors bundled on each end |
| X   | X   | X    | 5192963-2   | Run E2020 PGR-J22 to PEN-J17, 17000mm               | 2.4       | 61 | 0.82  | 21 | 7.0  | 178 | 6          | N/A | LFNC             | N/A                 | 85°             | optical fibers in conduit, 3 connectors bundled on each end |
| X   | X   | X    | 5194040     | Run E3023 SPW-J117 to MON-J8, 17000mm               | 2.5       | 64 | 0.44  | 11 | 6.0  | 152 | 25         | 22  | CL3, FT4         | 300                 | 80°             |                                                             |
| X   | X   | X    | 5194540     | Run P2026 OW to Mag-J11, Scan Room Display, 40000mm | 1.2       | 30 | 0.32  | 8  | 4.8  | 122 | 4          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 2 duplex connectors bundled on each end     |
| X   | X   | X    | 5194540-2   | Run P2026 OW to Mag-J11, Scan Room Display, 44000mm | 1.2       | 30 | 0.32  | 8  | 4.8  | 122 | 4          | N/A | OFNR             | N/A                 | 85°             | optical fibers, 2 duplex connectors bundled on each end     |
| X   | X   | X    | 5196200     | Run E4005 RF COM GND STUD to PGR-GND STUD, 10000mm  | 0.807     | 20 | 0.608 | 15 | 7.0  | 178 | 1          | 1/0 | UL AWM 1284, TEW | 600                 | 105°            | single conductor ground cable                               |
| X   | X   | X    | 5196200-2   | Run E4005 RF COM GND STUD to PGR-GND STUD, 17000mm  | 0.807     | 20 | 0.608 | 15 | 7.0  | 178 | 1          | 1/0 | UL AWM 1284, TEW | 600                 | 105°            | single conductor ground cable                               |

Illustration 6-10: System Cable List (Page 8)

|     |     |      | Part Number                          | Description                                             | Connector |       | Cable |    | Bend |     | Cond-uctors | AWG | Certifi-cation   | Volt Rating (volts) | Temp Rating (C) | Notes                         |
|-----|-----|------|--------------------------------------|---------------------------------------------------------|-----------|-------|-------|----|------|-----|-------------|-----|------------------|---------------------|-----------------|-------------------------------|
|     |     |      |                                      |                                                         | In        | mm    | In    | mm | In   | mm  |             |     |                  |                     |                 |                               |
| 450 | 750 | 450w |                                      |                                                         |           |       |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5195200-3                            | Run M4005 RF COM GND STUD to MAG-GND STUD, 13500mm      | 0.807     | 20    | 0.608 | 15 | 7.0  | 178 | 1           | 1/0 | UL AWM 1284, TEW | 600                 | 105°            | single conductor ground cable |
| X   | X   | X    | 5195200-4                            | Run M4005 RF COM GND STUD to MAG-GND STUD, 17500mm      | 0.807     | 20    | 0.608 | 15 | 7.0  | 178 | 1           | 1/0 | UL AWM 1284, TEW | 600                 | 105°            | single conductor ground cable |
| X   | X   | X    | 5195200-5                            | Run E4009 RF COM GND STUD to SPW-GND STUD, 2800mm       | 0.807     | 20    | 0.608 | 15 | 7.0  | 178 | 1           | 1/0 | UL AWM 1284, TEW | 600                 | 105°            | single conductor ground cable |
| X   | X   | X    | 5195200-6                            | Run E4010 RF COM GND STUD to PEN-GND STUD, 2800mm       | 0.807     | 20    | 0.608 | 15 | 7.0  | 178 | 1           | 1/0 | UL AWM 1284, TEW | 600                 | 105°            | single conductor ground cable |
| X   | X   | X    | 5196045                              | Run E4002, PEN-GND to PGR-GND, Ground Cable, 10000mm    | 0.531     | 13    | 0.18  | 5  | 5.0  | 127 | 1           | 10  | CMG              | 300                 | 80°             | single conductor ground cable |
| X   | X   | X    | 5196045-2                            | Run E4002, PEN-GND to PGR-GND, Ground Cable, 17000mm    | 0.531     | 13    | 0.18  | 5  | 5.0  | 127 | 1           | 10  | CMG              | 300                 | 80°             | single conductor ground cable |
| X   | X   | X    | 5196045-4                            | Run E4007, HEC-GND to PGR-GND, Ground Cable, 10000mm    | 0.531     | 13    | 0.18  | 5  | 5.0  | 127 | 1           | 10  | CMG              | 300                 | 80°             | single conductor ground cable |
| X   | X   | X    | 5196045-5                            | Run E4008, OW-GOC-GND to PGR-GND, Ground Cable, 30000mm | 0.531     | 13    | 0.18  | 5  | 5.0  | 127 | 1           | 10  | CMG              | 300                 | 80°             | single conductor ground cable |
| X   | X   | X    | 5197321-2                            | Run M3023 SPW-J117 to Mag-FJ2, 14000mm                  | 2.5       | 64    | 0.44  | 11 | 6.0  | 152 | 25          | 22  | CL3, FT4         | 300                 | 80°             | single conductor ground cable |
| 450 | 750 | 450w | Part Number                          | Description                                             | OD in.    | mm    |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5199506-2, -3, 46-282666P15, 16, 17, | 4" ID Air Hose                                          | 4.5       | 113.8 |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5213158                              | 3" ID Air Hose                                          | 3.4       | 86.3  |       |    |      |     |             |     |                  |                     |                 |                               |
|     |     | X    | 46-282666P18                         | 1.5" ID Air Hose                                        | 1.8       | 44.9  |       |    |      |     |             |     |                  |                     |                 |                               |
|     |     |      | 5264083,                             |                                                         |           |       |       |    |      |     |             |     |                  |                     |                 |                               |
|     |     |      | 5264086,                             |                                                         |           |       |       |    |      |     |             |     |                  |                     |                 |                               |
|     |     |      | 5264105,                             |                                                         |           |       |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5264033,                             | 1" Coolant Hose                                         | 1.3       | 32.6  |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5284039                              | 1/2" Coolant Hose                                       | 0.8       | 19.8  |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5225963                              | 4" Hose Insulation                                      | 5.5       | 139.2 |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5305915                              | 3" Hose Insulation                                      | 4.4       | 111.7 |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5225988                              | 1" Hose Insulation                                      | 2.3       | 58    |       |    |      |     |             |     |                  |                     |                 |                               |
| X   | X   | X    | 5197470                              | 1/2" Hose Insulation                                    | 2.3       | 57.9  |       |    |      |     |             |     |                  |                     |                 |                               |

## 4 Facility Supplied System Interconnects Specifications

The following table lists the required facility supplied system interconnects. Refer to [Illustration 6-11](#) for additional information.

**Table 6-6: Facility Supplied System Interconnects**

| Group | Between Units                         |                                                | Comments                                                                                                                                                                                                                     | Requirements                                                                    |
|-------|---------------------------------------|------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------|
|       | From                                  | To                                             |                                                                                                                                                                                                                              |                                                                                 |
| C01   | Facility Power                        | MDP                                            | Facility Power and Ground                                                                                                                                                                                                    | <a href="#">Chapter 2, MR Suite Electrical Requirements</a>                     |
| C02   | MDP                                   | HEC                                            | HEC Power                                                                                                                                                                                                                    |                                                                                 |
| C03   | MDP                                   | PGR                                            | PDU Power                                                                                                                                                                                                                    |                                                                                 |
|       | Facility Cooling Water                | HEC                                            | Cooling Water Supply                                                                                                                                                                                                         | <a href="#">Chapter 2, MR System Facility Water Requirements</a>                |
|       | Facility Cooling Water                | HEC                                            | Cooling Water Return                                                                                                                                                                                                         |                                                                                 |
| C04   | Facility Network                      | MON                                            | Facility must provide network/telephone access for the Magnet Monitor (MON) and Global Operator Cabinet (GOC). The MON connection must be available at all times.                                                            | <a href="#">Chapter 4, Magnet Monitor (MON) Requirements and Specifications</a> |
|       | Facility Network                      | GOC                                            |                                                                                                                                                                                                                              |                                                                                 |
| C05   | MDP                                   | E-Off Switch                                   | Facility must supply cable from MDP to E-Off Switch in Equipment Room                                                                                                                                                        | <a href="#">Chapter 2, MR Suite Electrical Requirements</a>                     |
|       | Facility Power                        | Outlet near MON                                | Facility outlet for MON power                                                                                                                                                                                                | <a href="#">Chapter 4, Magnet Monitor (MON) Requirements and Specifications</a> |
|       | Facility Power                        | MRU                                            | Facility power to MRU                                                                                                                                                                                                        | <a href="#">Chapter 3, Magnet Room Equipment Specifications</a>                 |
| C06   | GE Healthcare Supplied Cable from SPW | E-Off Switch in Control room or Equipment room | Facility must supply additional wiring between the GE Healthcare supplied cable and the E-Off switch if the length needed is greater than the usable length listed in <a href="#">MR System Interconnects Specifications</a> | <a href="#">Chapter 2, MR Suite Electrical Requirements</a>                     |

**Table 6-7: Optional Facility Supplied System Interconnects**


| Group | Between Units  |                 | Comments                                            | Requirements                                                                    |
|-------|----------------|-----------------|-----------------------------------------------------|---------------------------------------------------------------------------------|
|       | Facility Power | Outlet near MRE | Facility outlet for MRE power (optional, not shown) | <a href="#">Chapter 4, Magnetic Resonance Elastography (MRE) Specifications</a> |

## Illustration 6-11: Facility Supplied System Interconnects

## Facility Supplied System Interconnects

**Notes:**

1. Illustration is not to scale and component positioning/interconnect runs are typical
2. Solid lines are facility required interconnects; dashed lines are GE Healthcare supplied interconnects shown for clarity
3. Only GE Healthcare equipment interconnects are shown. Additional facility interconnects are required for non-GE Healthcare equipment (e.g., Magnet room DC Lighting)

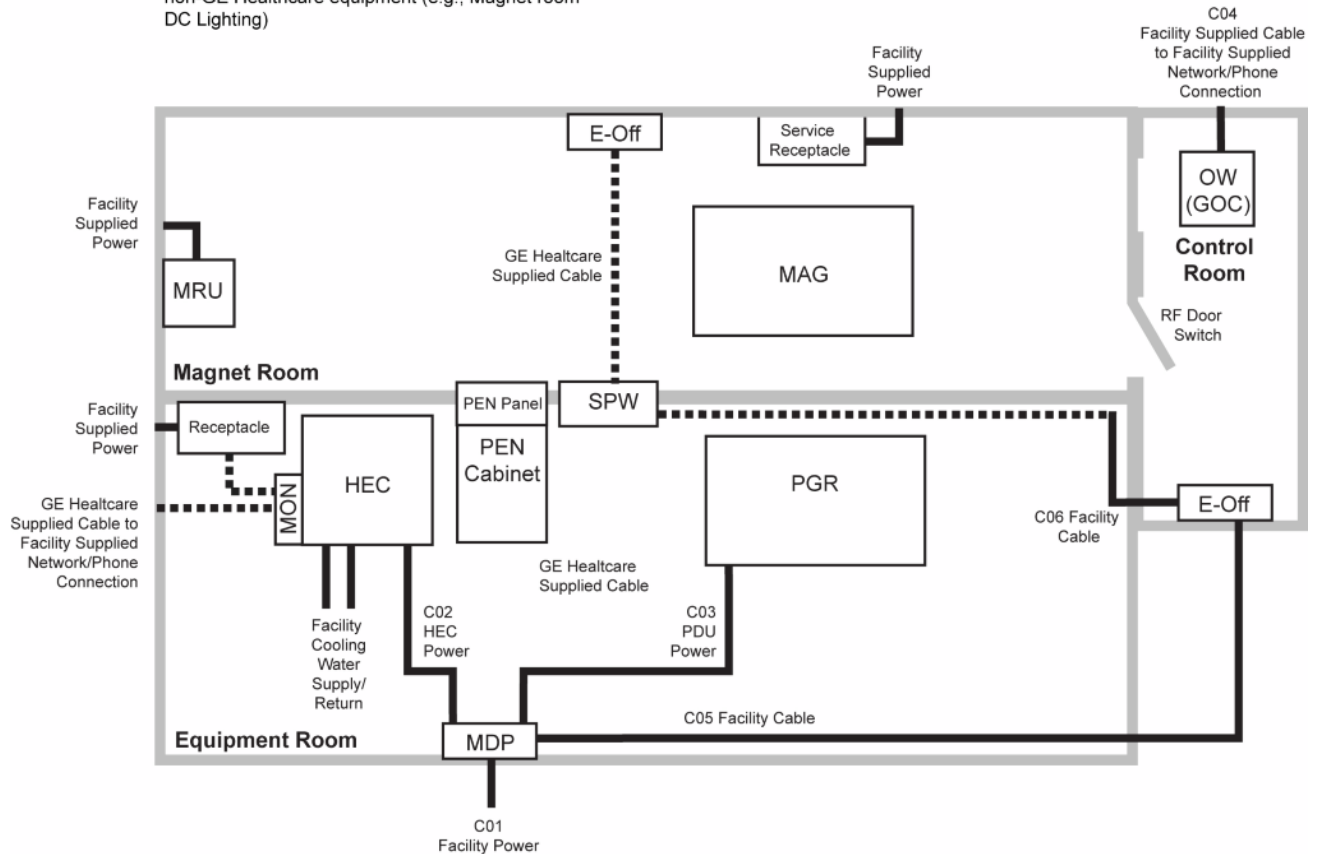

**NOTE:** GE Healthcare recommends installing the RF Door switch on the outside wall of the Magnet room.

# Chapter 7 Appendices

## 1 Glossary

### **BB**

Abbreviation for Broadband

### **Cryogen**

A substance for producing low temperatures. Liquid helium is the cryogen used to cool the magnet to approximately 4 Kelvin (-269°C or -452°F).

### **Cryostat**

An apparatus maintaining a very low constant temperature. The cryostat consists of one concentric, cylindrical container housed in an outer vacuum tight vessel. The magnet and shim coils are mounted in the inner container. The container is filled with liquid helium. The shields surrounding the inner container are kept cold by a refrigeration device.

### **Dewar**

A container with an evacuated space between two highly reflective walls used to keep low temperature substances at near-constant temperatures. Liquid helium is usually stored and shipped in dewars.

### **Exclusion Zone**

Area where the magnetic flux density is greater than five gauss. Personnel with cardiac pacemakers, neurostimulators and other biostimulation devices must NOT enter this zone. Signs are posted outside the five gauss line alerting personnel of this requirement. Since the magnetic field is three-dimensional, signs are also posted on floors above and below the Magnet Room in which the five gauss line exists.

### **Ferrous Material**

Any substance containing iron which is strongly attracted by a magnetic field.

### **Gauss (G)**

A unit of magnetic flux density. The earth's magnetic field strength is approximately one half gauss to one gauss depending on location. The internationally accepted unit is the tesla (1 Tesla = 10,000G and 1 milli Tesla = 10G).

### **Gradient**

The amount and direction of the rate of change in space of the magnetic field strength. In the magnetic resonance system, gradient amplifiers and coils are used to vary the magnetic field strength in the x, y, and z planes.

**Homogeneity**

Uniformity. The homogeneity of the static magnetic field is an important quality of the magnet.

**Isocenter**

Center of the imaging volume ideally located at the magnet center.

**Isogauss Line**

An imaginary line or a line on a field plot connecting identical magnetic field strength points.

**Magnetic Field (B)**

A condition in a region of space established by the presence of a magnet and characterized by the presence of a detectable magnetic force at every point in the region. A magnetic field exists in the space around a magnet (or current carrying conductor) and can produce a magnetizing force on a body within it.

**Magnetic Resonance (MR)**

The absorption or emission of electromagnetic energy by nuclei in a static magnetic field, after excitation by a suitable radio frequency field.

**Magnetic Shielding**

Using material (e.g. steel) to redistribute a magnetic field , usually to reduce fringe fields.

**NB**

Abbreviation for Narrow Band

**Quench**

Condition when a superconducting magnet becomes resistive thus rapidly boiling off liquid helium. The magnetic field reduces rapidly after a quench.

**Radio Frequency (RF)**

Frequency intermediate between audio frequency and infrared frequencies. Used in magnetic resonance systems to excite nuclei to resonance. Typical frequency range for magnetic resonance systems is 5-130 Mhz.

**Radio Frequency Shielding**

Using material (e.g. copper, aluminium, or steel) to reduce interference from external radio frequencies. A radio frequency shielded room usually encloses the entire magnet room.

**Resonance**

A large amplitude vibration caused by a relative small periodic stimulus of the same or nearly the same period as the natural vibration period of the system. In magnetic resonance imaging, the

radio frequency pulses are the periodic stimuli which are at the same vibration period as the hydrogen nuclei being imaged.

**Security Zone**

Area within the Magnet Room where the magnet is located. Signs are posted outside the Magnet Room warning personnel of the high magnetic field existing in the Magnet Room and the possibility of ferrous objects becoming dangerous projectiles within this zone.

**Shield Cooler Coldhead**

An external refrigeration device which maintains the shields inside the cryostat at a constant temperature.

**Shim Coils**

Shim coils are used to provide auxiliary magnetic fields in order to compensate for inhomogeneities in the main magnetic field due to imperfections in the manufacturing of the magnet or affects of steel in the surrounding environment.

**Shimming**

Correction of inhomogeneity of the main magnetic field due to imperfections in the magnet or to the presence of external ferromagnetic objects.

**Superconducting Magnet**

A magnet whose magnetic field originates from current flowing through a superconductor. Such a magnet is enclosed in a cryostat.

**Superconductor**

A substance whose electrical resistance essentially disappears at temperatures near zero Kelvin. A commonly used superconductor in magnetic resonance imaging system magnets is niobium-titanium embedded in a copper matrix.

**Tesla**

The internationally accepted unit of magnetic flux density. One tesla is equal to 10,000 gauss. One milli Tesla is equal to 10 gauss.

## 2 MR Site Vibration Test Guidelines

### 2.1 Test Measurements

1. Vibration measurements must be in the range of  $10^{-6}$  g. Test equipment must have the required sensitivity to these levels
2. Instrumentation must have a low tolerance to temperature effects since many times the low frequency thermal drift may influence the measurements
3. All measured data must be acquired real time. Recording of vibration data will not allow for a proper site survey, specifically when studying transient vibration and when searching for specific vibration sources
4. All analyses must be narrow-band Fast Fourier Transforms (FFT) over the frequency bands listed in [Table 7-1](#)
5. Time histories of the vibration must be recorded as acceleration levels vs. time. The resolution of the time history must be adjusted to clearly capture the transient event. The analyzer set-up will be site dependent and, in special cases, vibration response dependent. It is the responsibility of the vibration consultant to study the transient environment, capture data to confirm that transient activity exceeds the trigger level, then expand the time history data to exhibit the structural response

**Table 7-1: Frequency Bands for FFT**

| Frequency Band | Frequency Resolution  |
|----------------|-----------------------|
| 0.2 to 50 Hz   | $\Delta f = 0.125$ Hz |

### 2.2 Equipment (Spectral Analyzer) Set-Up

1. Frequency average should be a minimum of 20 linear averages (Do not use peak hold or 1/3 octave analysis)
2. Average and store should be a minimum of 10 plots to support the consistency of the site vibrations
3. Hanning windows must be applied to the entire spectra
4. Spectrum analyzers capable of these measurements are readily available for purchase or rental. Models, such as the HP 3560A, Nicolet Phaszer, B&K Pulse, and HP 35670, are all capable of making the site vibration measurements. Accelerometers must have the capability to measure from 0.2 Hz beyond 50 Hz. Time histories can be recorded using any of the analyzers listed above

**NOTE:** The equipment mentioned is for example only. It is the responsibility of the Engineering Test Firm to provide equipment that will allow measurements compliant with this guideline.

### 2.3 Data Collection

#### 2.3.1 Ambient Baseline Condition

1. All of the measurements listed above must be made in a “quiet” environment--i.e., areas where excessive traffic, subway trains, etc. do not exist. A vibration measurement must also be made



during periods without traffic or during periods of light traffic. Measurements must define the lowest levels of vibration possible at the site

2. The source of any steady state vibration, whose level exceeds the magnet specifications found in [Chapter 3, Magnet Room Structural Requirements](#), must be identified. A second measurement should be made with all of the identified contributors powered down if possible. In situations where it is not possible to power down equipment, vibration data must be collected to identify the specific source of the vibration concern. The majority of steady state vibration problems can be negated by isolating the vibration source

### 2.3.2 Normal Condition


1. All of the vibration measurements listed above must be repeated during periods of “normal” environmental conditions, including the FFTs and time histories. The transient measurements must be provided to define the dynamic disturbances the MR system may be exposed to. Transient analysis is required for a true assessment of the site
2. Special attention must be paid to the site assessment during the entire analysis. Since transient vibration is not easily addressed once the MR suite is fully constructed, the test consultant must fully understand the needs for this analysis. The source of any transient vibration must be identified and supported with vibration plots. If the source of any transient vibration is not locatable, it is recommended that the customer have an alternate location identified and the vibration studied
3. Transient vibration can be difficult to assess if the details are not understood. The **0.0005g, zero-to-peak trigger level** is a starting point to understanding the vibration stability. The transient vibration peak amplitude, structural (time variant) response, decay rate and an estimate of the number of events per unit of time would constitute a complete transient analysis. All transient failures must be supported by time history plots. The plots must clearly show the structural response, the frequency of the signature and the decay rate. From this data, GE Healthcare can help determine compliance with the vibration guidelines
4. The test consultant must provide site data to show the design recommendations for all sites/ building structures meet the Magnet specifications found in [Chapter 3, Magnet Room Structural Requirements](#)

## 2.4 Presentation/Interpretation of Results

1. The recommended format for site vibration data collection, presentation, and analysis is demonstrated in the examples in [Illustration 7-1](#), [Illustration 7-2](#), [Illustration 7-3](#), and [Illustration 7-4](#). Presentation of the data in any other format (linear units only) may result in incorrect interpretation and diagnosis of the site. Additional data collection or presentation methods are at the option of the vibration testing service
2. All plots must be properly annotated with:
  - a. Instrumentation setup including number of averages, frequency resolution, etc.
  - b. Test location
  - c. Test conditions:
    - i. Steady state
    - ii. Transient
    - iii. Heel drop

- iv. Normal environment
  - v. Typical traffic
  - vi. Any other conditions necessary to demonstrate understanding of potential sources of vibration
- 3. The customer's vibration testing service is responsible for interpreting the results and determine if that site meets GE Healthcare specifications
- 4. If the vibration levels are too high, additional data acquisition may be necessary to:
  - a. Determine the source of the vibration
  - b. Propose a solution to the problem
  - c. Find an alternate site location
- 5. Any questions regarding test equipment requirements, test parameters, or general questions should be discussed with the GE Healthcare Project Manager of Installation (PMI)

Illustration 7-1: Example of Site Environmental Vibration Levels

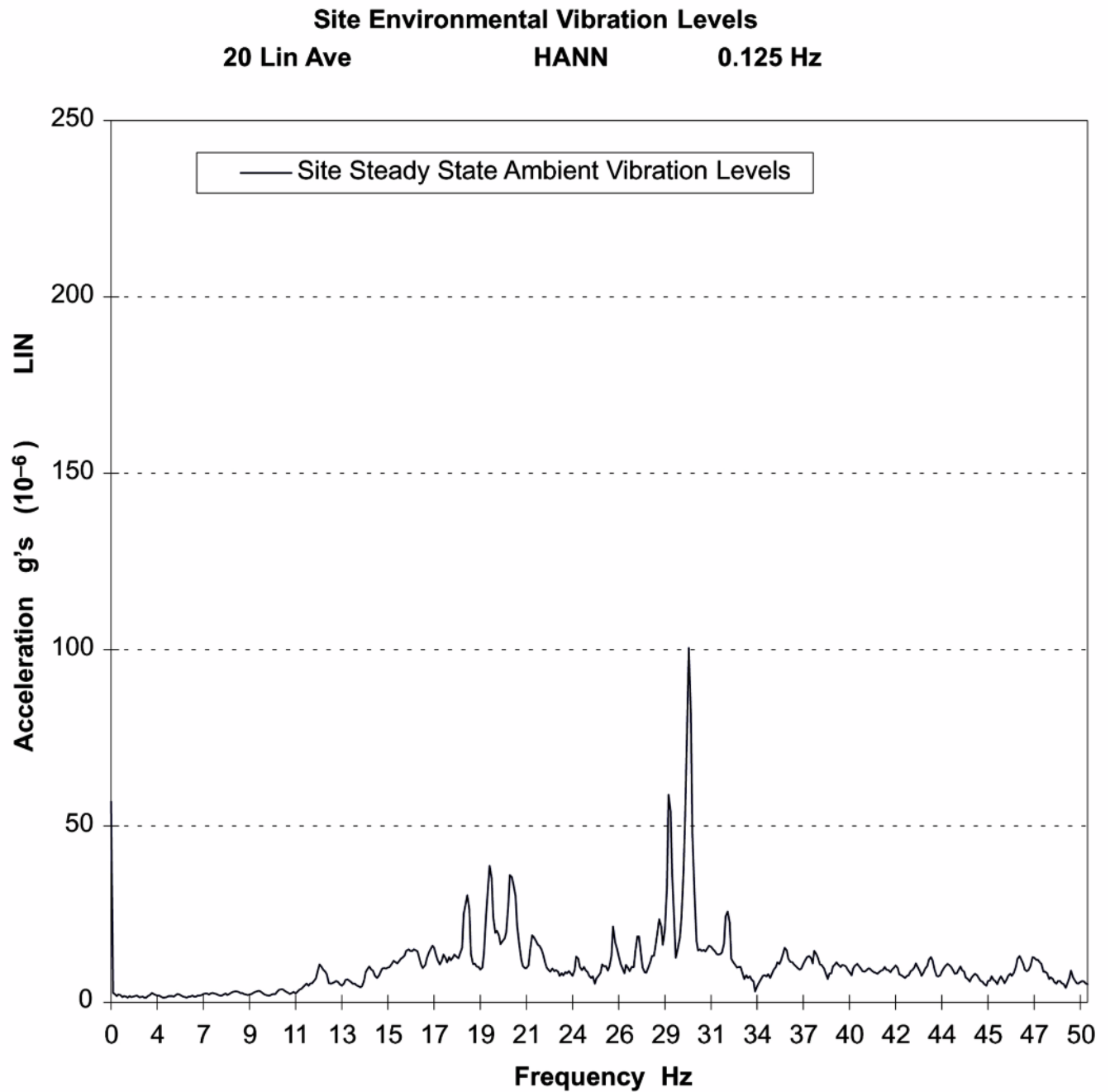

## Illustration 7-2: Example Site Environmental Vibration

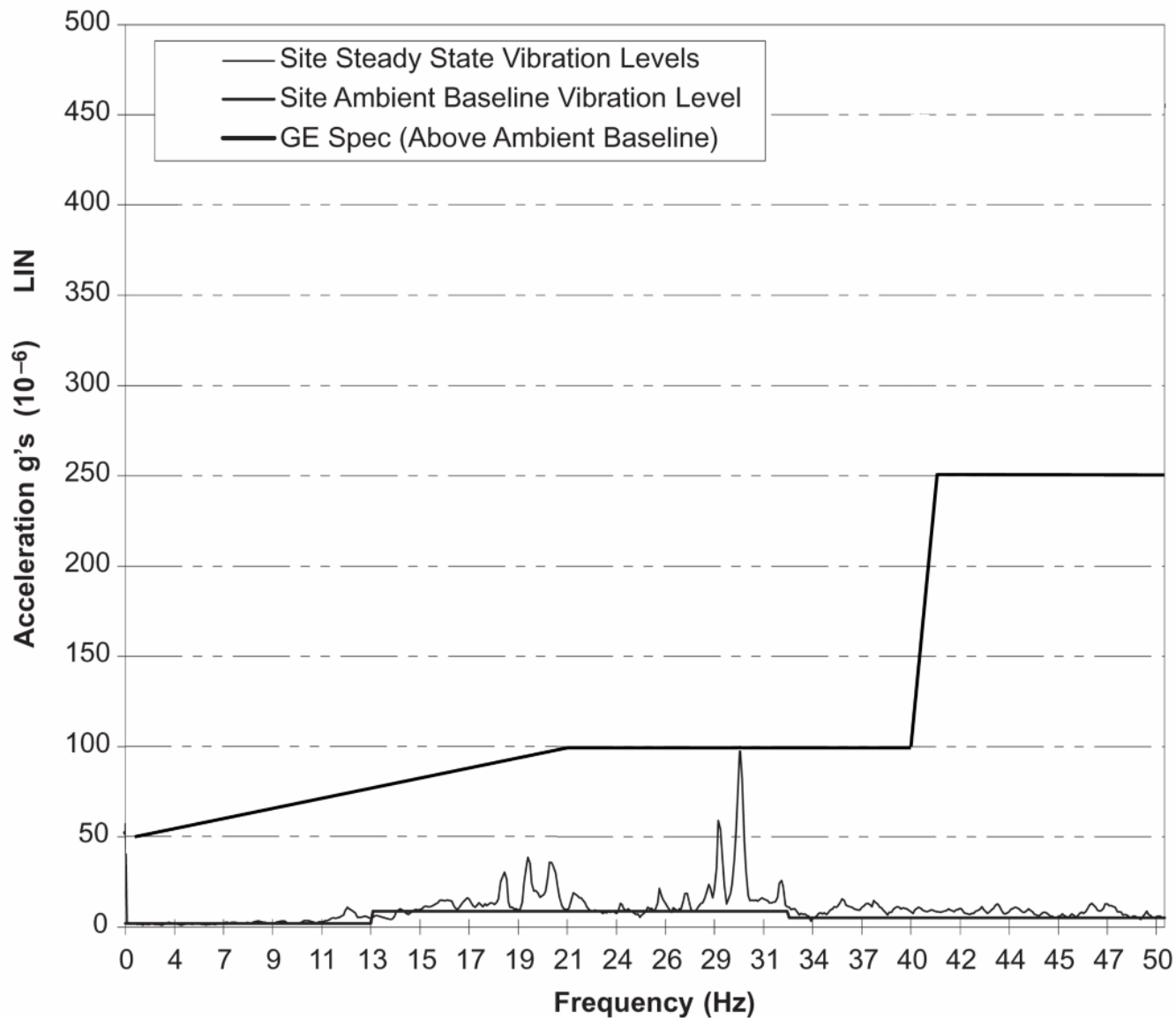
**EXAMPLE: Site Environmental Vibration vs. GE Spec. for 1.5T Magnet****20 Lin Ave****Hann****0.125Hz**

Illustration 7-3: Acceleration Time History


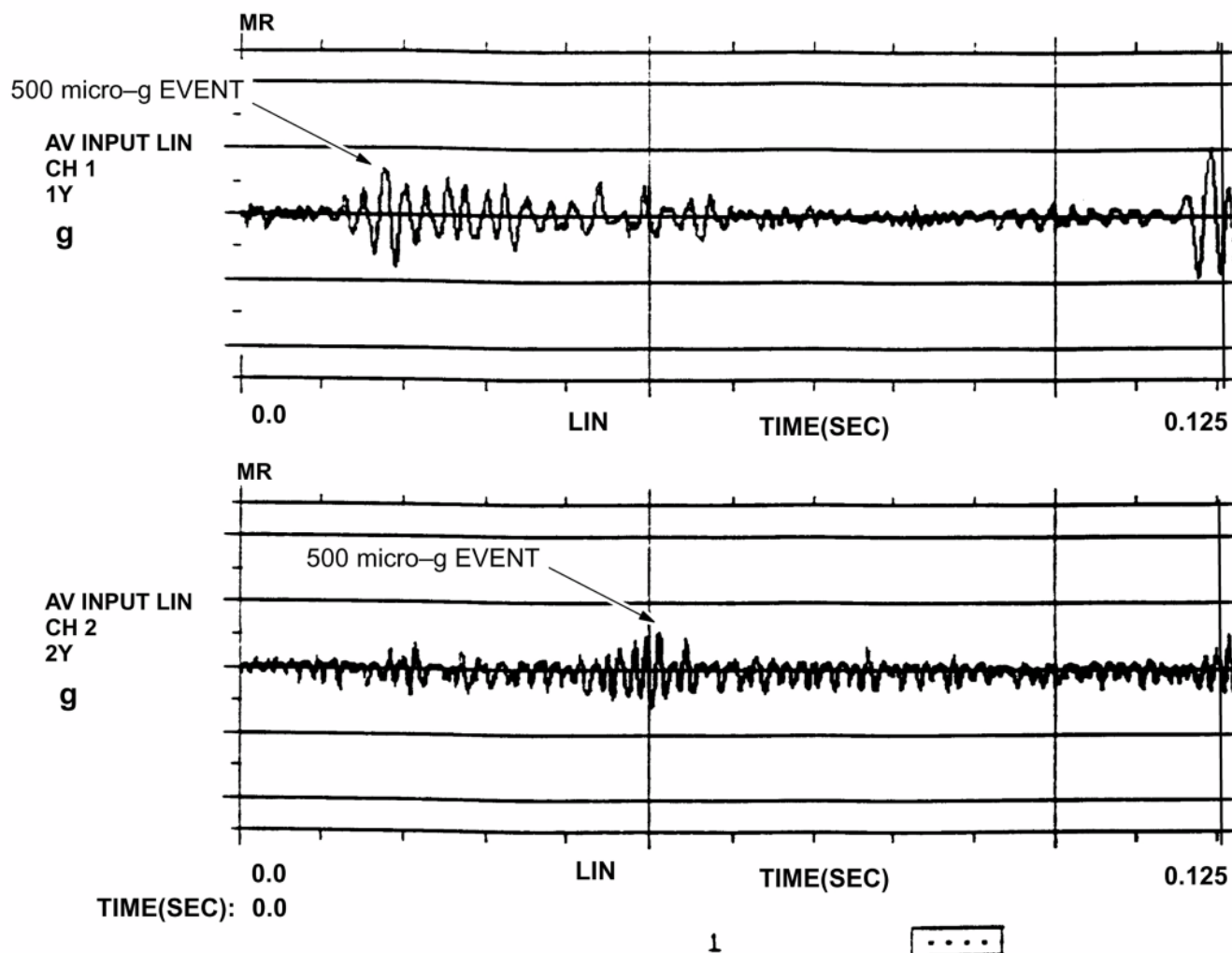


Illustration 7-4: Acceleration Time History (Zoomed In on Transient Event)

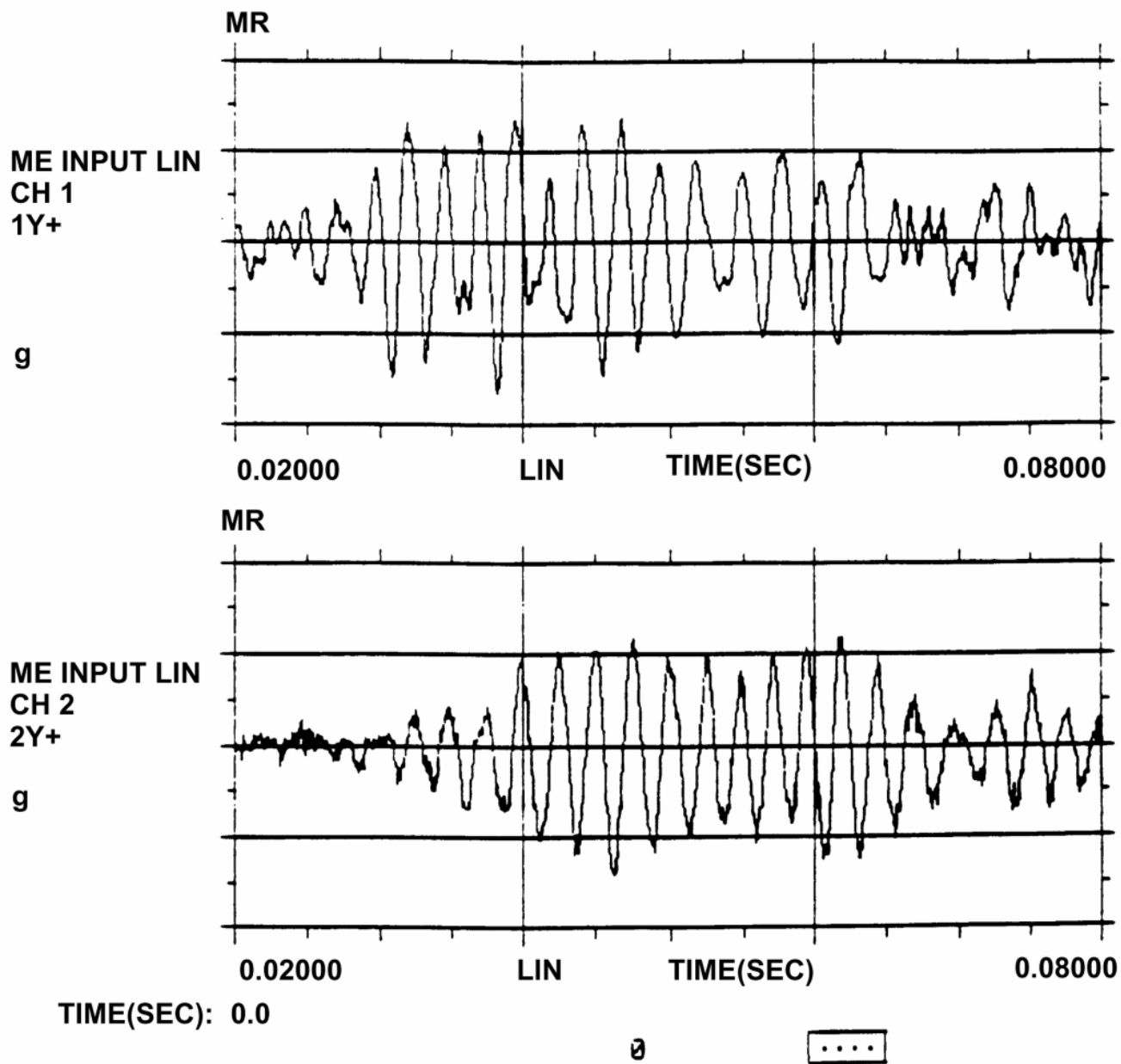

## 3 RF Shielding Background and Guidelines

### 3.1 RF Shielding Background

RF sources which can adversely affect image quality may be generated by discrete frequency or broadband noise (RF) sources.

#### 3.1.1 Discrete Frequency

Discrete RF interferences are narrowband and are fixed frequency in nature. The Magnet Room must be RF shielded from RF sources so external RF energy does not degrade the MR system RF receivers at the system imaging frequencies. Some potential sources for discrete frequency signals are radio station transmitters, mobile or hand-held RF transmitting devices, etc.

#### 3.1.2 Broadband RF Noise

Broadband RF noise is a single transient or continuous series of transient disturbances caused by an electrical discharge. Low humidity environmental conditions will have higher probability of electrical discharge. The electrical discharge can occur due to electrical arcing (micro arcing) or merely a static discharge. Some potential sources capable of producing electrical discharge include:

- Loose hardware/fasteners vibration or movement (electrical continuity must always be maintained)
- Flooring material including raised access flooring (panels & support hardware) and carpeting
- Electrical fixtures
  - Lighting fixtures
  - Track lighting
  - Emergency lighting
  - Battery chargers
  - Outlets
- Ducting for HVAC and cable routing
- RF Shield seals (walls, doors, windows, etc.).

### 3.2 RF Shielded Room Grounding

- Primary Ground: All RF Shield components (walls, floor, ceiling, etc) must be electrically bonded together to form one common ground plane which is connected to the Facility Grounding Conductor.
- Secondary Ground (Other grounds that connect the outside of the RF Shield Room to earth grounds are called secondary grounds): Secondary grounds must not compromise Ground Isolation of 1,000 ohms or greater.

## 4 RF Shielded Enclosure Test Guidelines

### 4.1 Introduction

This document describes the procedure and methodology of performing an RF shielding effectiveness verification testing on enclosures that will house GE Healthcare Magnetic Resonance Imaging (MRI) equipment. MRI equipment is sensitive to RF energy from sources outside of the shielded enclosure. To ensure proper operation of the MRI equipment, the shielded enclosure must attenuate local RF signals to levels that do not cause interference.

**NOTE:** RF Shielding Performance is based on plane-wave measurements. *H* and *F* field tests are not required, but are allowed as needed for diagnostic purposes.

#### 4.1.1 Purpose of Test Plan

The purpose of this test plan is to describe a series of RF shielding effectiveness tests to demonstrate compliance of an MRI shielded enclosure to the requirements of GE Healthcare.

The test procedure described in this guideline is a modification of MIL-STD-285 and IEEE Std 299-1991. This procedure provides a thorough evaluation of the shield integrity at the upper end of the frequency range of interest showing any RF leakage that may cause imaging problems. These testing guidelines ensure that the electromagnetic environment inside of the enclosure will meet the requirements of GE Healthcare.

#### 4.1.2 Not used

### 4.2 Applicable Documents

|                   |                                                                                                                                              |
|-------------------|----------------------------------------------------------------------------------------------------------------------------------------------|
| MIL-STD-285       | MILITARY STANDARD ATTENUATION MEASUREMENTS FOR ENCLOSURES, ELECTRO-MAGNETIC SHIELDING, FOR ELECTRONIC TEST PURPOSES, METHOD OF; 25 June 1956 |
| IEEE Std 299-1991 | IEEE STANDARD FOR MEASURING THE EFFECTIVENESS OF ELECTROMAGNETIC SHIELDING ENCLOSURES; 2 July 1991                                           |

### 4.3 Test Sample Set-Up

The shielded enclosure under testing will be set up in a normal configuration that consists of:

- Magnet installed including all floor mounting bolts
- RF shielded door(s)
- Waveguide penetrations, HVAC, vents, medical gas lines, etc.
- AC power supplied through low-pass filters
- Patient view window, skylights, windows, hatches, etc.
- Blank penetration panel installed, dimensionally equivalent to the GE panel and the same mounting hardware to be used with the GE penetration panel.

For safety reasons, the enclosure will be electrically grounded during the shielding effectiveness test. Any variances from the normal configuration will be noted in the certification report.

### 4.4 Shielding Effectiveness

This test procedure determines the worst case shielding effectiveness based on the lowest test point reading obtained. The lowest reading obtained will be the reading of the room.



## 4.5 Measurement Procedure

To simulate the effects of external RF sources, the transmitting antenna will be located outside the enclosure on a plane parallel to the face of the enclosure wall at a distance of 6 ft. (1.8 m) unless physically constrained to a lesser separation. The areas of least effectiveness are located by searching the inside of the enclosure with the antenna connected to the spectrum analyzer.

### 4.5.1 Test Position

The transmitting antenna will be positioned in front of all critical areas (doors, windows, filters, penetration areas, etc.) at a minimum of every 20 ft. (6.1 m) of the wall. The receiving antenna is scanned over all panel section joints (where accessible) at the floor, wall, and ceiling for a minimum of 10 ft. (3.05 m) in all directions from the location of transmitting antenna. The receiving antenna will be at a minimum of 1 ft. (0.3 m) from the shield. For areas that are inaccessible for direct location of the transmitting antenna, the inside of that area will still be scanned using the receiving antenna with the transmitting antenna positioned in front of the adjacent wall or test position.

### 4.5.2 Frequency Range

The standard frequency for shielding measurements will be 100 MHz  $\pm$  10 MHz (150 MHz  $\pm$  10 MHz for 3T). This allows the frequency to be adjusted slightly to avoid interference from local active transmitters and/or RF noise from other sources. Test frequency utilized will be noted in the certification report.

### 4.5.3 Free Field Calibration

The incident field (free field) is measured by the following procedure:

Position the transmitting antenna parallel to the exterior wall of the enclosure at a distance of 6 ft. (1.8 m) using horizontal polarization, unless physically constrained to a lesser separation, in which case a separate reference will be established and documented at the new test distance. The receiving antenna will be placed between the transmitting antenna and 1 ft. (0.3 m) from the exterior wall of the enclosure. The receiving antenna will be moved vertically and horizontally to achieve maximum signal strength. The receiving antenna will be placed no closer than 2 in. (51 mm) from the exterior wall of the enclosure and in line with the transmitting antenna. The maximum received voltage at the test frequency will be recorded.

## 4.6 Enclosure Power Reference Isolation



### NOTICE

This section does not apply to upgrades.

To prevent personal hazard, it is necessary for the enclosure to be properly grounded.

To minimize common mode currents, the ungrounded enclosure should be isolated from the ground with a minimum of 1000 ohms of DC resistance. The isolation measurement is performed by the following procedure:

- All power to the enclosure is removed. For safety reasons, an AC voltage measurement will be made to verify that no power is connected.

- With electrical power and intentional ground disconnected, connect the test instrument between the shielded enclosure and AC power ground.
- Take a reading and record the value.

**NOTE:** This test must be made using either an isolated, current limited, high-voltage (>150 VDC) DC source and DMM to read drop across the limiting resistor or a Megger instrument capable of reading values less than 1000 ohms. Conventional resistance meters employing test sources of 9 VDC or less will not be used.

## 4.7 Test Equipment

Test equipment will be selected to provide measuring capabilities as described in these testing guidelines. The signal source, amplifier, antennas, and receiver or spectrum analyzer will be such that the difference between the induced reference voltage and the receiver sensitivity is at least 6 dB greater than the required attenuation specification.

The signal source and power amplifier will output a CW signal for a nominal test frequency of 100 mHz (150 mHz for 3T). The receiver or spectrum analyzer and preamplifier (if required) will provide adequate sensitivity to permit attenuation measurements to be made at the specified limits. Dipole antennas and other miscellaneous equipment required to transmit and receive the proper RF fields will be used.

The absolute performance calibration of the equipment requiring calibration will be performed on an as-needed basis in accordance with MIL-STD-45662. The calibration period will not exceed one year. The test equipment tolerances of at least  $\pm 2\%$  frequency and  $\pm 2$  dB amplitude will be met. Equipment certifications will be traceable to the National Institute of Standards and Technology (NIST). All equipment will be verified for proper operation between and after each series of tests by repeating the reference readings at the specified frequency(s).

## 4.8 Data Recording and Verification

Measurements will be performed by qualified responsible EMC test personnel. The test must be performed in the presence of a GE Healthcare representative unless other arrangements were made by GE Healthcare. All data collected during the course of the tests will be recorded on standardized data sheets. The data sheets will include the test location, frequency, reference level, measured enclosure level, and attenuation level.

## 4.9 Test Report

A final certification report will be provided after the test is performed. This report will include all recorded data necessary for the evaluation of the shielded enclosure test results and will list any changes pertinent to the test set-up or shielding effectiveness. The certification report will also include the test procedures and a list of the actual equipment used during the test. Along with the data sheet, there will be a presentable drawing showing the shape of the enclosure, all test point locations, doors, filters, windows, and existing building walls.

## 5 Acoustic Background and Design Guidelines

The acoustic information is provided for site planning and architectural design activities to address acoustics to meet local regulations and customer requirements. For more information about recommended safety procedures regarding patient exposure to MR generated acoustic levels, see the MR Safety Guide included with the system Users Manual.

### 5.1 Acoustic Background

A typical MR suite has two types of acoustic noise issues. The first is the acoustics within the rooms in which the patients and technicians are impacted by the noise of the MR system as the gradients are pulsed. The second is noise transmitted to other spaces via airborne and structureborne paths.

#### 5.1.1 Airborne

The airborne transmission path entails the excitation of air within the magnet room; the resonator module consisting of the magnet, RF coil, and gradient coil generates acoustic noise similar to an intense loud speaker. The airborne noise passes through walls via any openings, i.e. small holes, cracks, HVAC ducts, and waveguides, into surrounding spaces within and possibly beyond the confinements of the building. Acoustic energy can transmit across distances of significant length.

Examples of airborne acoustics issues may include the following (not limited to only these) :

- MR Operator exposure at Operator Workstation (i.e. Operator viewing in-line with the patient inside the magnet may require a higher acoustic attenuation window)
- Image reading rooms adjacent to Magnet Room, may be separated by hallways
- Secretarial, offices, meeting rooms, patient rooms (ICU, exam, primary care, etc.)
- Adjacent residential areas/spaces
- In-house library facilities

#### 5.1.2 Structureborne

The structureborne transmission path is the result of mechanical excitation of the floor/building structure causing the building to vibrate. The vibration of the surfaces at surrounding spaces then radiates as acoustic noise. Acoustic energy can transmit across distances of significant length.

**NOTE:** Less than 5% of installed base sites have experienced structureborne acoustic issues.

Examples of structureborne acoustics issues may include the following (not limited to only these):

- Areas directly above or below the Magnet Room, may not always be an issue
- Image reading rooms adjacent to Magnet Room, may be separated by hallways
- Secretarial, offices, meeting rooms, patient rooms (ICU, exam, primary care, etc.)
- Adjacent residential areas/spaces
- In-house library facilities

### 5.2 Acoustic Design Guidelines

#### 5.2.1 Magnet Room

Noise generated by the MR system is inherent to the operation of the system. The sound quality (human perception) within the Magnet Room can be modified by including sound absorbing materials to make the room sound more subdued and less harsh. The measured sound levels via

a sound level meter will not change. However, the measured sound levels can be reduced only when the sound level generated by the MR System is reduced.

Sound quality improvements can be achieved by the following:

- Use ceiling tiles with fiberglass panels having a 2 inch (51mm) thickness set into the standard T-bar grid system.
- Adding fiberglass panels to the side walls covering approximately 20% of the side wall surface area. The panels should focus on covering the top half of the side walls. Panels could take many different and decorative shapes to improve the sterile look of the rooms. Typically panels might be on the order of 4ft x 6ft (1.2m x 1.8m) with a thickness of 4 inches (102mm) or equivalent. Panels shape could vary to produce mosaic effects to meet the customer preference. Any decorative materials used to cover the wall panels must be porous so that sound waves can pass through with ease. In principle, a person should be able to breath through the material with ease. Fire retardant cloth should be used. The NRC (Noise Reduction Coefficient) of the panels should be 0.95 or better when mounted against a hard surface such as drywall or concrete.

### 5.2.2 Inter-Spacial Areas

Acoustic Noise Control to mitigate noise from being transmitted to other spaces often amounts to paying attention to small details while working with ordinary construction materials. The key objectives are to eliminate all cracks and gaps in the wall construction while making sure that the doors, walls, floor, and ceiling have adequate transmission loss via mass or special double wall construction along with good fitting massive doors.

The entire Magnet must be surrounded by walls with substantial mass and/or double wall construction so that noise is contained in the room and not allowed to pass through into nearby spaces. Wall junctions must be sealed with acoustical sealant so that noise waves do not escape from the room. In principle, if the room were filled with smoke and under a positive pressure, no smoke would leak from the room.

#### 5.2.2.1 Wall Construction

Wall Construction will entail ordinary building materials in a careful configuration.

- The preferred wall would have an ASTM STC 50 construction which entails the use of standard wall construction of steel studs (typically 3-5/8 inch (92 mm)) with 2 layers of Type X drywall (typically 5/8 (16 mm)) on each side totaling 4 layers and fiberglass batt in the stud cavity. All drywall must be overlapped by 6 inches (152 mm) or more. Beads of (USG) Acoustical Caulking (non-hardening) would be used around the entire perimeter of the drywall. Any form of wall penetration should be avoided. Any necessary wall penetrations must be sealed using combination of Acoustical Caulking (non-hardening) and fiberglass batt material. See examples of wall construction shown in [Illustration 7-5](#) and [Illustration 7-6](#).
- The top of the wall must join the ceiling/floor above so that no cracks or gaps occur. If metal pan is used on the ceiling/floor (above), then flute seals would be used to seal the gaps between the drywall and the pan. Alternately drywall can be cut out to fit into the flutes. Acoustical caulking (non-hardening) will be used to seal the remaining cracks and gaps.

Illustration 7-5: Example Of Wall Construction For Airborne Noise Control - Option 1


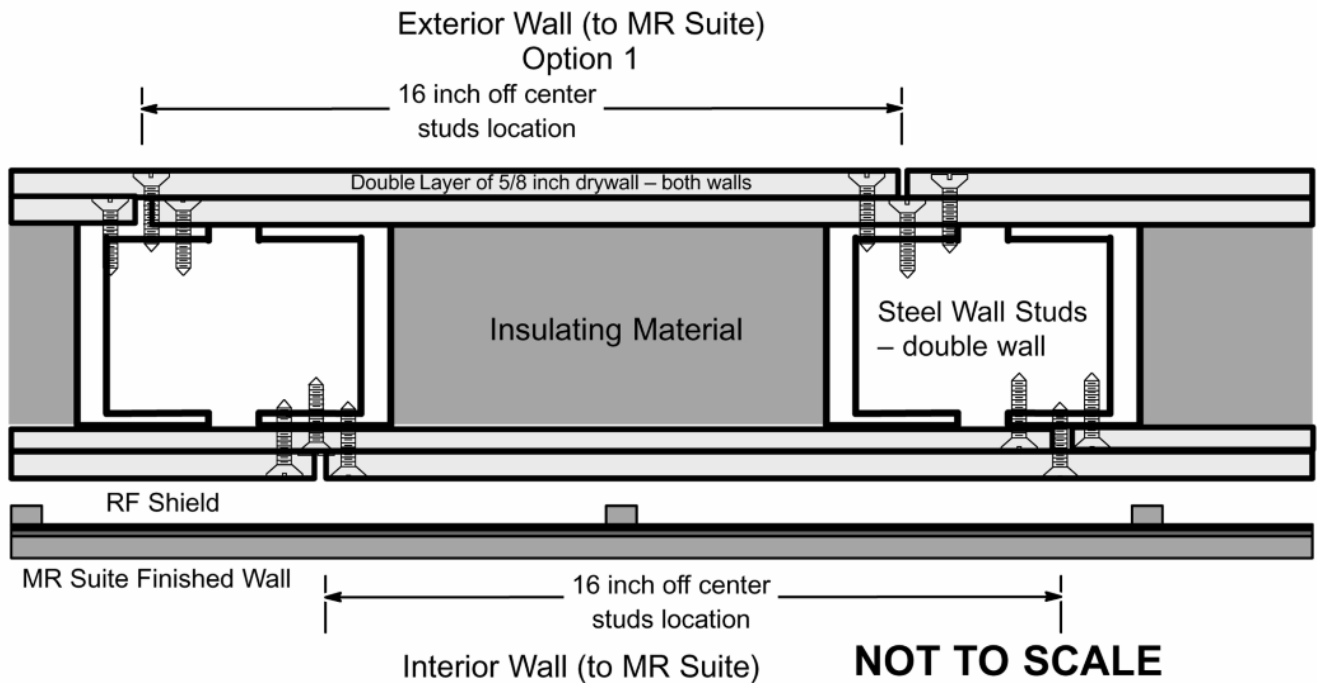


Illustration 7-6: Example Of Wall Construction For Airborne Noise Control - Option 2

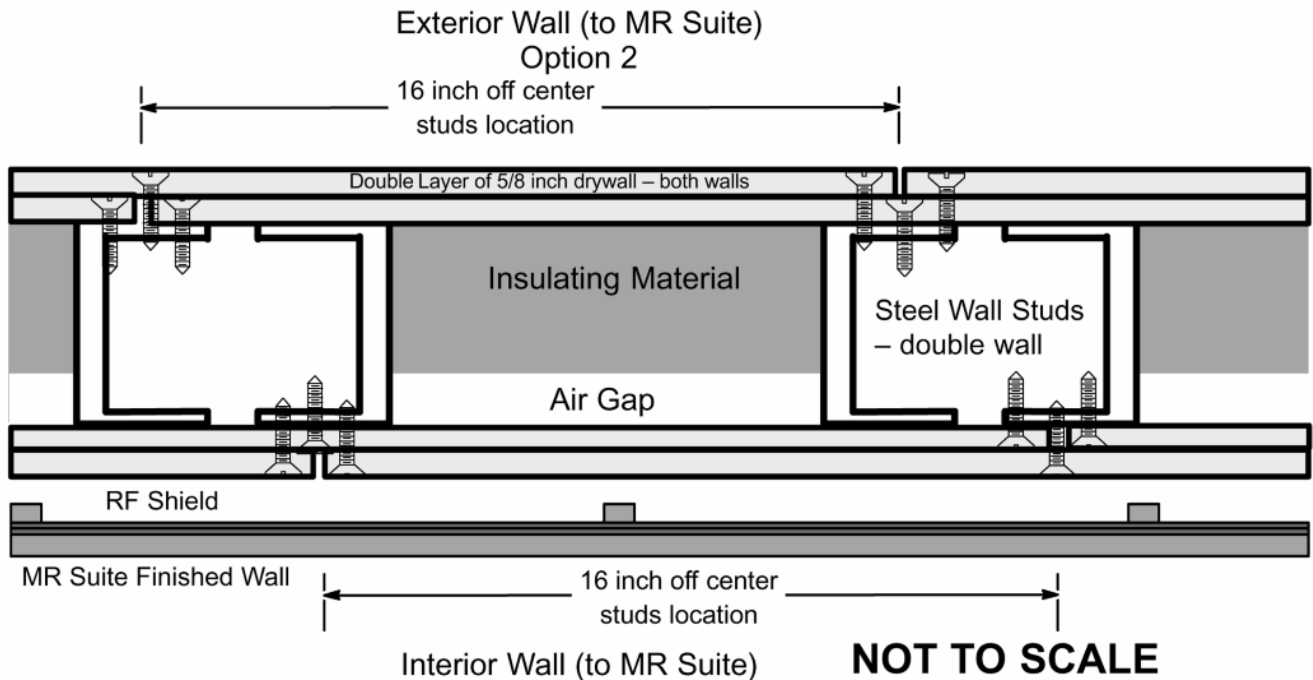

### 5.2.2.2 High Bay RF Room

A high bay RF Room is a self contained RF Room which has open air space between the RF Room ceiling and the building floor above. The air space is an acoustic transmission path. Acoustic energy must be reduced to minimize this transmission of energy through this path.

In cases where the Magnet is to be installed in a high bay, it may be most effective to enclose the RF Room with its own drywall and steel stud room. The key difference being a ceiling assembly that mimics the sidewall construction to contain noise.

- Normal high STC stud walls from above would be used to support a ceiling assembly constructed of structural C channel with two layers of drywall on each side (total of 4 layers) with fiberglass batt in the cavity.
- Penetrations should be avoided via the use of surface mounted lights. HVAC and ducts passing through the ceiling, party wall or side walls would require acoustic noise attenuation in the form of inline silencers. Gaps and cracks would be sealed between the ceiling, party wall or vertical side walls and the cryogen vent plumbing. In essence the Magnet would be enclosed in a drywall "doghouse".

#### 5.2.2.3 Miscellaneous Plumbing, RF Windows and RF Doors

Other construction details are equally important to mitigate noise transmission to meet the intended goal.

- Pipes (gas or water) and electrical conduit or Magnet Room signal cables must be sealed where they penetrate the walls or ceiling. A heavy mastic material such as Duxseal™ is appropriate.
- RF windows should be purchased as window/frame units with an STC rating obtained from laboratory testing per ASTM standards. STC 50 to 60 windows are needed. The installation must include proper sealing to avoid sound leaks.
- RF doors should be selected to provide an STC 50 to 60 to quell the noise. Contact RF Shield Room supplier for selection of RF doors that meet the local acoustic codes and site acoustic requirements. RF door seals must be selected to prevent small gaps around the door perimeter and at the door threshold. RF door seals would either require periodic replacement or a door seal that would last the life of the Magnet Room.

## 6 Sample Calculation AC Power Equipment Minimum Distance

This is a sample calculation to determine minimum distance from a feeder, transformer, or other AC electrical source, using the Formula found in [Chapter 2, MR Suite Magnetic Field Specifications](#) Electrical Current subsection to determine minimum distance from a feeder, transformer, or other AC electrical source.

$$I \text{ (amps)} = 20X^2 \text{ (meters)} \div S \text{ (meters)}$$

Note that the formula has 3 variables, if you have 2 of them, you can calculate the 3rd. In this example, we calculate the minimum distance  $X$  from the source, in this case a main electrical feeder carrying 450 amps of current in a 5 inch conduit.


$$x^2 \text{ (meters)} = \{I \text{ (amps)} \times S \text{ (Meters)}\} / 20$$

$$x = \sqrt{\frac{\{I \text{ (amps)} \times S \text{ (Meters)}\}}{20}}$$

Rearranging:

Note that the separation "S" is the spacing between the conductors and when all 3 conductors are run in a single conduit, "S" is simply the diameter of the conduit.

$$S = 5 \text{ inches} = 0.127 \text{ meters}$$

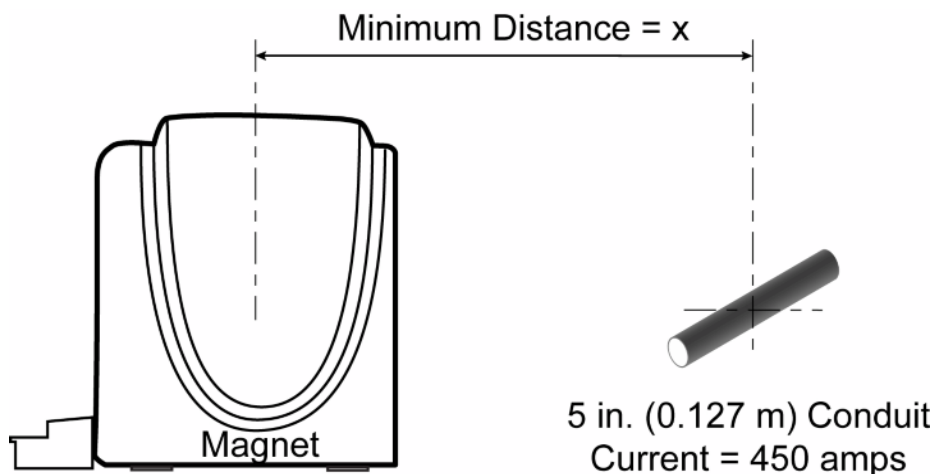

$$x = \sqrt{\frac{\{450 \text{ (amps)} \times 0.127 \text{ (Meters)}\}}{20}}$$

The conduit should be 1.7 meters or 5.6 feet from the magnet isocenter.

In other situations, the spacing "S" may be the spacing between HV feeders, the distance between transformer lugs, or the spacing between conduits when the phase conductors are run in separate conduits.

What if it is too close? Keep in mind that if this is an existing condition, you should request an *EMI study* to quantify the magnitude and direction of the AC disturbances. The calculation is worse

case and does not take into account the vector direction of the AC interference. The magnet is only sensitive to AC disturbances that are directed horizontally (magnet z-axis). Also the calculation does not account for any magnetic shielding effect of steel conduit.



## 7 Selecting Magnet Anchor Size

The following is an example to illustrate the selection of proper anchors to install a Magnet in a building with 2000 psi (13.8 MPa) concrete. For this example the area is not under seismic requirements.

1. Determine magnet clamping force (for the Magnet: 2500 lbs + 200 lbs = 2700 lbs or 11,100 N + 900 N = 12,000 N).
2. Refer to the examples of anchor vendor catalogs below to select anchor diameter and embedment which meets the clamping force (tension) determined in Step 1.  
Diameter : Min. 0.625 inch Max. 1.25 inch  
For 8 inch embedment select 3/4 inch diameter  
For 4.5 inch embedment select 1 inch diameter  
or  
Diameter : Min. M16 Max. M32  
For 130 mm embedment select M20 diameter  
For 114 mm embedment select M24 diameter
3. The vendor instructions and torque to the maximum recommended level for the anchor selected in Step 2 must be provided to the RF Shield Room vendor for proper installation of the anchor and equipment.

**Illustration 7-7: Allowable Anchor Loads in Concrete (English Units)**

| ANCHOR<br>DIAMETER<br>in. (mm)<br>See Note 1                                                                     | EMBEDMENT<br>DEPTH<br>in. (mm) | 2000 psi (13.8 MPa) |                  | 3000 psi (20.7 MPa) |                  | 4000 psi (27.6 MPa) |                  | 6000 psi (41.4 MPa) |                  |
|------------------------------------------------------------------------------------------------------------------|--------------------------------|---------------------|------------------|---------------------|------------------|---------------------|------------------|---------------------|------------------|
|                                                                                                                  |                                | TENSION<br>lb (kN)  | SHEAR<br>lb (kN) | TENSION<br>lb (kN)  | SHEAR<br>lb (kN) | TENSION<br>lb (kN)  | SHEAR<br>lb (kN) | TENSION<br>lb (kN)  | SHEAR<br>lb (kN) |
| 5/8<br>(15.9)                                                                                                    | 2 3/4<br>(70)                  | 1250<br>(5.6)       | 2800<br>(12.5)   | 1600<br>(7.1)       | 3070<br>(13.7)   | 1810<br>(8.1)       | 3330<br>(14.8)   | 1920<br>(8.5)       | 3330<br>(12.5)   |
|                                                                                                                  | 4<br>(102)                     | 1870<br>(8.3)       | 3330<br>(14.8)   | 2400<br>(10.7)      | 3330<br>(14.8)   | 2930<br>(13.0)      | 3330<br>(14.8)   | 3200<br>(14.2)      | 3330<br>(12.5)   |
|                                                                                                                  | 7<br>(178)                     | 2500<br>(11.2)      | 3330<br>(14.8)   | 3010<br>(13.4)      | 3330<br>(14.8)   | 3650<br>(16.2)      | 3330<br>(14.8)   | 3650<br>(16.2)      | 3330<br>(12.5)   |
| 3/4<br>(19.1)                                                                                                    | 3 1/4<br>(83)                  | 1550<br>(6.9)       | 2880<br>(12.8)   | 1950<br>(8.7)       | 3310<br>(14.7)   | 2350<br>(10.5)      | 3730<br>(16.6)   | 2610<br>(11.6)      | 4800<br>(21.4)   |
|                                                                                                                  | 4 3/4<br>(121)                 | 2510<br>(11.2)      | 4510<br>(20.1)   | 3250<br>(14.5)      | 4650<br>(20.7)   | 3870<br>(17.2)      | 4800<br>(21.4)   | 4670<br>(20.8)      | 4800<br>(21.4)   |
|                                                                                                                  | 8<br>(203)                     | 2930<br>(13.0)      | 4800<br>(21.4)   | 3870<br>(17.2)      | 4800<br>(21.4)   | 4530<br>(20.2)      | 4800<br>(21.4)   | 5120<br>(22.8)      | 4800<br>(21.4)   |
| 1<br>(25.4)                                                                                                      | 4 1/2<br>(114)                 | 3120<br>(13.9)      | 6080<br>(27.0)   | 3870<br>(17.2)      | 6770<br>(30.1)   | 4610<br>(20.5)      | 7470<br>(33.2)   | 4800<br>(21.4)      | 7470<br>(33.2)   |
|                                                                                                                  | 6<br>(152)                     | 4400<br>(19.6)      | 7470<br>(33.2)   | 6400<br>(28.5)      | 7470<br>(33.2)   | 7200<br>(32.0)      | 7470<br>(33.2)   | 7330<br>(32.6)      | 7470<br>(33.2)   |
|                                                                                                                  | 9<br>(229)                     | 5600<br>(24.9)      | 7470<br>(33.2)   | 8000<br>(35.59)     | 7470<br>(33.2)   | 9390<br>(41.77)     | 7470<br>(33.2)   | 9390<br>(41.8)      | 7470<br>(33.2)   |
| <b>Note 1</b> All shaded values fail to meet the clamping force (tension), therefore are not acceptable anchors. |                                |                     |                  |                     |                  |                     |                  |                     |                  |

Illustration 7-8: Allowable Anchor Loads in Concrete (Metric Units)

| ANCHOR<br>DIAMETER<br>See Note 1                                                                                 | EMBEDMENT<br>DEPTH<br>mm (in.) | 13.8 MPa (2000 psi) |                  | 20.7 MPa (3000 psi) |                  | 27.6 MPa (4000 psi) |                  | 41.4 MPa (6000 psi) |                  |
|------------------------------------------------------------------------------------------------------------------|--------------------------------|---------------------|------------------|---------------------|------------------|---------------------|------------------|---------------------|------------------|
|                                                                                                                  |                                | TENSION<br>kN (lb)  | SHEAR<br>kN (lb) | TENSION<br>kN (lb)  | SHEAR<br>kN (lb) | TENSION<br>kN (lb)  | SHEAR<br>kN (lb) | TENSION<br>kN (lb)  | SHEAR<br>kN (lb) |
| M16                                                                                                              | 105<br>(4 1/8)                 | 11.2<br>(2500)      | 25.1<br>(5650)   | 20.9<br>(4705)      | 39.9<br>(8965)   | 24.2<br>(5450)      | 10125<br>(45.0)  | 6900<br>(30.7)      | 10550<br>(46.9)  |
| M20                                                                                                              | 130<br>(5 1/8)                 | 25.1<br>(5650)      | 52.9<br>(11900)  | 30.7<br>(6910)      | 58.7<br>(13195)  | 36.4<br>(8175)      | 14490<br>(64.5)  | 10005<br>(44.5)     | 14490<br>(64.5)  |
| M24                                                                                                              | 155<br>(6 1/8)                 | 30.0<br>(6735)      | 61.2<br>(13760)  | 36.9<br>(8300)      | 70.5<br>(15855)  | 43.9<br>(9860)      | 29.8<br>(17950)  | 57.7<br>(12980)     | 95.6<br>(21490)  |
| <b>Note 1</b> All shaded values fail to meet the clamping force (tension), therefore are not acceptable anchors. |                                |                     |                  |                     |                  |                     |                  |                     |                  |

## 8 Magnet Cryogenic Venting Pressure Drop Reference Tables

Use the following tables to calculate the cryogenic vent pressure drop through the pipe used.

**Table 7-2: 1.5T Magnet Cryogenic Vent System Pressure Drop Matrix**

| Inside dia. of vent pipe (D) | Distance of vent system component from magnet |            | Pressure drop for straight pipe with smooth inside surface |       | Pressure Drop Per Elbow Used Anywhere Within A 20 Ft (6.1 M) Vent Segment |       |                                |       |                                  |       |                                  |       |
|------------------------------|-----------------------------------------------|------------|------------------------------------------------------------|-------|---------------------------------------------------------------------------|-------|--------------------------------|-------|----------------------------------|-------|----------------------------------|-------|
|                              |                                               |            |                                                            |       | Standard Sweep                                                            |       |                                |       | Long Sweep                       |       |                                  |       |
|                              |                                               |            |                                                            |       | Std sweep 45° elbow (K = 15 Ft)                                           |       | Std sweep 90° elbow (K = 30Ft) |       | Long sweep 45° elbow (K = 7.5Ft) |       | Long sweep 90° elbow (K = 15 Ft) |       |
|                              | ft                                            | m          | psi/ft                                                     | kPa/m | psi                                                                       | kPa   | psi                            | kPa   | psi                              | kPa   | psi                              | kPa   |
| 8 in. (203.2 mm)             | 00- 20                                        | 0 -6.1     | 0.10                                                       | 2.26  | 1.10                                                                      | 7.58  | 2.06                           | 14.20 | 0.55                             | 3.79  | 1.03                             | 7.10  |
|                              | 20- 40                                        | 6.1 -12.2  | 0.21                                                       | 4.75  | 2.10                                                                      | 14.48 | 3.70                           | 25.51 | 1.03                             | 7.10  | 1.85                             | 12.76 |
|                              | 40- 60                                        | 12.2 -18.3 | 0.30                                                       | 6.79  | 2.88                                                                      | 19.86 | 5.21                           | 35.92 | 1.44                             | 9.93  | 2.60                             | 17.92 |
|                              | 60-80                                         | 18.3-24.4  | 0.38                                                       | 8.60  | 3.70                                                                      | 25.51 | 6.71                           | 46.27 | 1.85                             | 12.76 | 3.36                             | 23.17 |
|                              | 80-100                                        | 24.4-30.5  | 0.47                                                       | 10.63 | 4.52                                                                      | 31.17 | 8.22                           | 56.68 | 2.26                             | 15.58 | 4.11                             | 28.34 |
| 10 in. (254.0 mm)            | 0- 20                                         | 0 - 6.1    | 0.03                                                       | 0.68  | 0.55                                                                      | 3.79  | 0.82                           | 5.65  | 0.27                             | 1.86  | 0.41                             | 2.83  |
|                              | 20- 40                                        | 6.1 -12.2  | 0.07                                                       | 1.58  | 0.82                                                                      | 5.65  | 1.51                           | 10.41 | 0.41                             | 2.83  | 0.75                             | 5.17  |
|                              | 40- 60                                        | 12.2 -18.3 | 0.10                                                       | 2.26  | 1.23                                                                      | 8.48  | 2.19                           | 15.10 | 0.62                             | 4.27  | 1.10                             | 7.58  |
|                              | 60- 80                                        | 18.3 -24.4 | 0.12                                                       | 2.71  | 1.51                                                                      | 10.41 | 2.74                           | 18.89 | 0.75                             | 5.17  | 1.37                             | 9.45  |
|                              | 80- 100                                       | 24.4 -30.5 | 0.16                                                       | 3.62  | 1.92                                                                      | 13.24 | 3.43                           | 23.65 | 0.96                             | 6.62  | 1.71                             | 11.79 |
| 12 in. (304.8 mm)            | 0- 20                                         | 0 - 6.1    | 0.013                                                      | 0.29  | 0.27                                                                      | 1.86  | 0.41                           | 2.83  | 0.14                             | 0.97  | 0.21                             | 1.45  |
|                              | 20- 40                                        | 6.1 -12.2  | 0.027                                                      | 0.61  | 0.41                                                                      | 2.83  | 0.82                           | 5.65  | 0.21                             | 1.45  | 0.41                             | 2.83  |
|                              | 40- 60                                        | 12.2 -18.3 | 0.041                                                      | 0.93  | 0.55                                                                      | 3.79  | 1.10                           | 7.58  | 0.27                             | 1.86  | 0.55                             | 3.79  |
|                              | 60- 80                                        | 18.3 -24.4 | 0.054                                                      | 1.22  | 0.69                                                                      | 4.76  | 1.37                           | 9.45  | 0.34                             | 2.34  | 0.69                             | 4.76  |
|                              | 80- 100                                       | 24.4 -30.5 | 0.069                                                      | 1.56  | 0.96                                                                      | 6.62  | 1.51                           | 10.41 | 0.48                             | 3.31  | 0.75                             | 5.17  |
|                              | 100-120                                       | 30.5-36.6  | 0.08                                                       | 1.81  | 1.09                                                                      | 7.52  | 1.77                           | 12.20 | 0.55                             | 3.79  | 0.88                             | 6.07  |
|                              | 120-140                                       | 36.6-42.7  | 0.10                                                       | 2.26  | 1.27                                                                      | 8.76  | 2.07                           | 14.30 | 0.63                             | 4.34  | 1.04                             | 7.17  |
|                              | 140-160                                       | 42.7-48.8  | 0.11                                                       | 2.49  | 1.43                                                                      | 9.86  | 2.36                           | 16.30 | 0.72                             | 4.96  | 1.19                             | 8.21  |
|                              | 160-180                                       | 48.8-54.9  | 0.12                                                       | 2.71  | 1.60                                                                      | 11.00 | 2.53                           | 17.40 | 0.80                             | 5.52  | 1.27                             | 8.76  |
|                              | 180-200                                       | 54.9-61.0  | 0.17                                                       | 3.85  | 1.75                                                                      | 12.10 | 2.93                           | 20.20 | 0.88                             | 6.07  | 1.47                             | 10.14 |
| 14 in. (355.6 mm)            | 0- 20                                         | 0 - 6.1    | 0.008                                                      | 0.055 | 0.20                                                                      | 1.38  | 0.301                          | 2.08  | 0.102                            | 0.70  | 0.15                             | 1.03  |
|                              | 20- 40                                        | 6.1 -12.2  | 0.017                                                      | 0.12  | 0.30                                                                      | 2.07  | 0.602                          | 4.15  | 0.154                            | 1.06  | 0.30                             | 2.07  |
|                              | 40- 60                                        | 12.2 -18.3 | 0.026                                                      | 0.18  | 0.40                                                                      | 2.76  | 0.808                          | 5.57  | 0.198                            | 1.37  | 0.40                             | 2.76  |
|                              | 60- 80                                        | 18.3 -24.4 | 0.034                                                      | 0.23  | 0.51                                                                      | 3.52  | 1.01                           | 6.96  | 0.250                            | 1.72  | 0.51                             | 3.52  |
|                              | 80- 100                                       | 24.4 -30.5 | 0.043                                                      | 0.30  | 0.71                                                                      | 4.90  | 1.11                           | 7.65  | 0.353                            | 2.43  | 0.55                             | 3.79  |
|                              | 100-120                                       | 30.5-36.6  | 0.050                                                      | 0.34  | 0.80                                                                      | 5.52  | 1.30                           | 8.96  | 0.40                             | 2.76  | 0.64                             | 4.41  |
|                              | 120-140                                       | 36.6-42.7  | 0.063                                                      | 0.43  | 0.933                                                                     | 6.43  | 1.52                           | 10.48 | 0.46                             | 3.17  | 0.76                             | 5.24  |
|                              | 140-160                                       | 42.7-48.8  | 0.069                                                      | 0.48  | 1.05                                                                      | 7.24  | 1.73                           | 11.93 | 0.52                             | 3.59  | 0.87                             | 6.00  |
|                              | 160-180                                       | 48.8-54.9  | 0.076                                                      | 0.52  | 1.18                                                                      | 8.14  | 1.85                           | 12.76 | 0.59                             | 4.07  | 0.93                             | 6.41  |
|                              | 180-200                                       | 54.9-61.0  | 0.11                                                       | 0.76  | 1.29                                                                      | 8.89  | 2.15                           | 14.82 | 0.64                             | 4.41  | 1.08                             | 7.45  |
| 16 in. (406.4 mm)            | 0- 20                                         | 0 - 6.1    | 0.0053                                                     | 0.037 | 0.153                                                                     | 1.05  | 0.230                          | 1.59  | 0.078                            | 0.54  | 0.115                            | 0.79  |
|                              | 20- 40                                        | 6.1 -12.2  | 0.013                                                      | 0.09  | 0.229                                                                     | 1.58  | 0.460                          | 3.17  | 0.118                            | 0.81  | 0.229                            | 1.58  |
|                              | 40- 60                                        | 12.2- 18.3 | 0.020                                                      | 0.14  | 0.306                                                                     | 2.11  | 0.618                          | 4.26  | 0.152                            | 1.05  | 0.306                            | 2.11  |
|                              | 60- 80                                        | 18.3 -24.4 | 0.026                                                      | 0.18  | 0.390                                                                     | 2.69  | 0.773                          | 5.33  | 0.191                            | 1.32  | 0.390                            | 2.69  |

| Inside dia. of vent pipe (D) | Distance of vent system component from magnet |            | Pressure drop for straight pipe with smooth inside surface |       | Pressure Drop Per Elbow Used Anywhere Within A 20 Ft (6.1 M) Vent Segment |      |                                             |       |                                               |      |                                               |      |
|------------------------------|-----------------------------------------------|------------|------------------------------------------------------------|-------|---------------------------------------------------------------------------|------|---------------------------------------------|-------|-----------------------------------------------|------|-----------------------------------------------|------|
|                              |                                               |            |                                                            |       | Standard Sweep                                                            |      |                                             |       | Long Sweep                                    |      |                                               |      |
|                              |                                               |            |                                                            |       | Std sweep 45° elbow (K = 15 F <sub>t</sub> )                              |      | Std sweep 90° elbow (K = 30F <sub>t</sub> ) |       | Long sweep 45° elbow (K = 7.5F <sub>t</sub> ) |      | Long sweep 90° elbow (K = 15 F <sub>t</sub> ) |      |
|                              | ft                                            | m          | psi/ft                                                     | kPa/m | psi                                                                       | kPa  | psi                                         | kPa   | psi                                           | kPa  | psi                                           | kPa  |
|                              | 80- 100                                       | 24.4 -30.5 | 0.033                                                      | 0.23  | 0.543                                                                     | 3.74 | 0.850                                       | 5.86  | 0.270                                         | 1.86 | 0.421                                         | 2.90 |
|                              | 100-120                                       | 30.5-36.6  | 0.038                                                      | 0.26  | 0.613                                                                     | 4.23 | 0.995                                       | 6.86  | 0.310                                         | 2.14 | 0.490                                         | 3.38 |
|                              | 120-140                                       | 36.6-42.7  | 0.048                                                      | 0.33  | 0.714                                                                     | 4.92 | 1.16                                        | 8.00  | 0.352                                         | 2.43 | 0.581                                         | 4.01 |
|                              | 140-160                                       | 42.7-48.8  | 0.052                                                      | 0.36  | 0.803                                                                     | 5.54 | 1.32                                        | 9.10  | 0.398                                         | 2.74 | 0.666                                         | 4.59 |
|                              | 160-180                                       | 48.8-54.9  | 0.058                                                      | 0.40  | 0.903                                                                     | 6.23 | 1.42                                        | 9.79  | 0.451                                         | 3.11 | 0.712                                         | 4.91 |
|                              | 180-200                                       | 54.9-61.0  | 0.084                                                      | 0.56  | 0.987                                                                     | 6.81 | 1.64                                        | 11.31 | 0.490                                         | 3.38 | 0.826                                         | 5.70 |

**Notes**

1. Elbows with angles greater than 90° **must not** be used.
2. The table data is based on the followings:
  - a. Initial flow conditions at magnet interface
  - b. Gas temperature starting at 4.5 Kelvin (-452° F or -268°C).
  - c. Helium gas flow rate of 2,737 cubic feet per minute (77.5 cubic meters per minute)
  - d. 45° standard sweep elbow K = 15 F<sub>t</sub>
  - e. 90° standard sweep elbow K = 30 F<sub>t</sub>
  - f. 45° long sweep elbow K = 7.5 F<sub>t</sub>
  - g. 90° long sweep elbow K = 15 F<sub>t</sub>



© 2009 General Electric Company.

GE Medical Systems, a General Electric Company, going to market as GE Healthcare.

3200 N. Grandview Boulevard

Waukesha, Wisconsin 53188

USA

[www.gehealthcare.com](http://www.gehealthcare.com)

