

Effectivity

The following CT Equipment is affected:

2276776 LightSpeed Plus Table Assembly

Purpose

The purpose of this FMI is to add an IV Pole to the CT table.

Related FMIs

- FMI 25280: Table Handle Installation
- FMI 25282: Table Cradle Installation

Furnished Materials

FMI Kit Part Number 2299386 includes the items listed below.

Qty	Description	Part #
1	FMI Document	2299188-100
1	IV Pole	2309994

People Power

Allow 0.10 hours of labor for the completion of this FMI.

Procedure

- 1.) Remove IV Pole from its packaging.
- 2.) Place assembled IV Pole in either the left or right accessory holes in the table cradle handle.

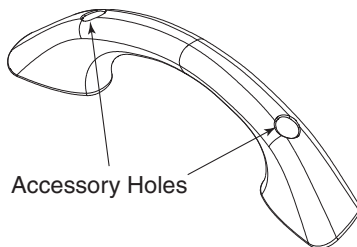


Figure 1 Accessory holes in table handle

Figure 2 (below) shows the IV pole installed on the LightSpeed Plus CT table.

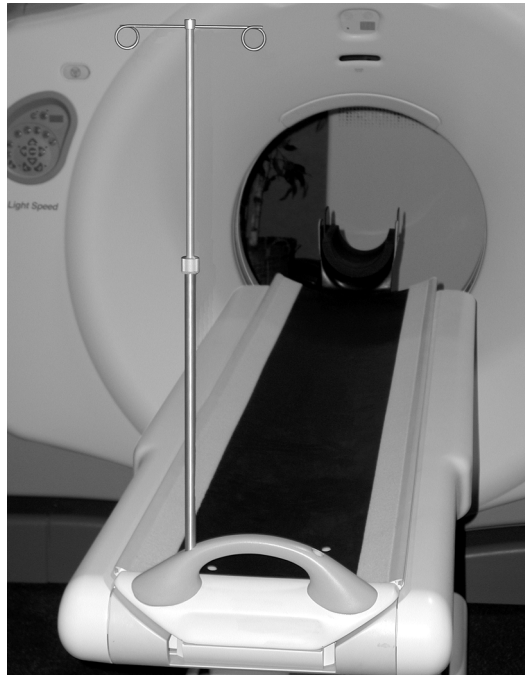


Figure 2 CT table, with IV Pole installed

End of Document