

**ADMINISTRATIVE INFORMATION****FMI No: 25311A****FMI Title: LS2002 M3 MDAS Software FMI**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS-AM:	Issue Date:
☐ Mandatory	GEMS-E:	Close Date:
☒ Mandatory on Request	GEMS-A:	

☐ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☐	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2353934-100****Do NOT leave unsecured on site.**



Completion Certification Report

Mail To

GEMS - Americas

GE Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe

GEMS - E / ESC
FMI Admin. 322
BP34
F 78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25311A FMI Title: LS2002 M3 MDAS Software FMI

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code:

- 1.) Modification Complete or previously modified.
- 2.) FMI not required per FMI instruction.
- 3.) Not modified. Unit in storage. FMI kit returned to stock.
- 4.) Modification refused or equipment altered. Customer signature and title required.
- 5.) Unit location unknown or unit not in area. Product Locator notified.

**Customer Refusal or Equipment Altered
- Modification Pending.**

Installer
Comments

Effectivity

1 "First Wave"

For the "first wave" of implementation, LightSpeed Ultra and Plus (MDAS) systems with the following Gantry system IDs are affected:

System ID	Gantry Model #	Gantry Serial Number
206292CT	2281177	00000282675CN8
206731MAINCT	2281177	00000284734CN1
208743SCTU	2281177	00000279658CN9
231935CT2	2281177-2	00000288650CN5
262785CT3	2281177	00000289220CN6
305996CT1	—	—
330996CTultra	2281177	00000291171CN7
336716BCT20	2281177	00000277304CN2
405271ULT	2281177	00000289403CN8
408358CT	2281177	00000273736CN9
415353CT4	2281177	00000289224CN8
415923PCT	2281177	00000276421CN5
502560ULER	2281177	00000289405CN3
502560ULT2	2281177	00000289225CN5
610645LSU	2281177	00000288983CN0
707525SLSU	2281177	00000276254CN0
713790MCT6	2281177	00000286459CN3
732240CT1	2281177	00000281829CN2
734712ct2	2281177	00000288098CN7
775885ULTRA	2281177	00000289679CN3
785623CT	2281177	00000277221CN8
805569CTLS1	2281177	00000278429CN6
856751ULTRA	2281177	00000281752CN6
910212092	—	—
910212095	—	—
925779CT	2281177	00000282140CN3
QY0009	2281177	00000276868CN7
QZ0001	2281177	00000276406CN6
QZ0006	2281177	00000277914CN8
QZ0010	2281177	00000277700CN1
QZ0013	2281177	00000278283CN7
QZ0022	2281177	00000281694CN0
YG0437	2281177	00000276423CN1

2 "Second Wave"

For the "second wave" of implementation, LightSpeed Ultra and Plus (MDAS) systems with the following Gantry system IDs are affected:

System ID	Gantry Model #	Gantry Serial Number
—	2281177	282936CN4
—	2281177-2	286554CN1
—	2281177-2	286628CN3
—	2281177-2	286870CN1
00069CTS01	2281177	00000284320CN9
01224CTS02	2281177	00000278606CN9
082421040072	—	—
082421040078	2281177	00000285393CN5
082421160059	2281177	00000273249CN3
190052CT02	—	—
260010CT08	2281177	00000284665CN7
260010CT09	2281177	00000285478CN4
315464CTLS2	2281177-2	289995CN3
417269LSU	2281177	291172CN5
475290CT04	2281177	00000279429CN5
502100CT06	2281177	00000279004CN6
606439LS	2281177	291504CN9
661949CT	2281177	277128CN5
734712CT1	2281177-2	287072CN3
786090CT03	2281177	00000288822CN0
805473CT	2281177-4	873727YM2
853210122	2281177-2	00000290203CN9
853210124	2281177	00000290649CN3
910212092	2281177	00000284063CN5
910216028	2281177-2	00000287472CN5
920219024	2281177	00000288918CN6
985867CTLS	2281177	289217CN2
A5105144	2281177	00000278570CN7
A5183747	2281177	00000278748CN9
B4192932	2281177	00000284985CN9
DPM99914	2281177	00000285496CN6
E4191718	2281177	00000286460CN1
HAUG02CT01	2281177	00000284735CN7
HC0074CT03	2281177	00000273800CN3
M2863012	2281177	00000281200CN6
M4016550	2281177	00000279104CN4
M40345140	2281177	00000279214CN1

System ID	Gantry Model #	Gantry Serial Number
M4169024	2281177-2	00000287452CN7
M9483895	2281177	00000281667CN6
M9886822	2281177	00000284429CN8
M9957416	2281177	00000288954CN1
MC6946CT05	—	—
NPV13916	2281177	00000273889CN6
P02002CT05	2281177	00000284777CN0
PCWU02CT02	2281177	00000278607CN7
PM1633CT03	2281177	00000278747CN1
QX0008	2281177	00000279182CN0
QY0015	2281177	00000279548CN2
QZ0005	2281177	00000269939CN5
QZ0007	2281177	00000278572CN3
QZ0008	2281177	00000281201CN4
QZ0009	2281177	00000279048CN3
QZ0011	2281177	00000278474CN2
QZ0015	2281177	00000281199CN0
QZ0017	2281177	00000281415CN0
QZ0020	2281177	00000281198CN2
QZ0023	2281177	00000285041CN0
QZ0025	2281177	00000282147CN8
QZ0028	2281177	00000281666CN8
UC1220CT06	2281177	00000277786CN0

Purpose

The purpose of this FMI is to support customers who are level 4 escalations and customers who were promised Cardia Burst and Smart Step.

Furnished Materials

FMI 25311 Kit (Part Number 2354014) includes the following items:

Qty	Description	Part #
1	FMI Document	2353934-100
1	LS2002 SW Set (MDAS) (see details, below)	2350061

The LS2002 SW Set (MDAS) (P/N 2350061) includes the following items:

Qty	Description	Part #
1	LS2002 H3 M3 OS and Apps (MDAS) (see details, below)	2348415
1	LS2002 Software Installation Manual	2348269-100
1	Application Tips and Workarounds	2296447-100
1	CT LightSpeed Ultra CD-ROM	2304119-246
1	Serial #, Model #	46-302200P7
1	440 mA Scan Option	2345411

The LS2002 MDAS M3 Software Kit (P/N 2348415) includes the following items:

Qty	Description	Part #
1	Operating System Disk 1	2305087-2
1	Operating System Disk 2	2305088
1	Operating System Disk 3	2305089
1	Applications Software	2348414
1	LS2002 M3 Patch CD	2349518


Required Tools

Qty	Description	Part #
1	LightSpeed Ultra/Plus (MDAS) Class C	2344424
1	LightSpeed Ultra/Plus (MDAS) Class M	2344425

People Power

4 hours labor and 1 hour of travel.

Steering Guide

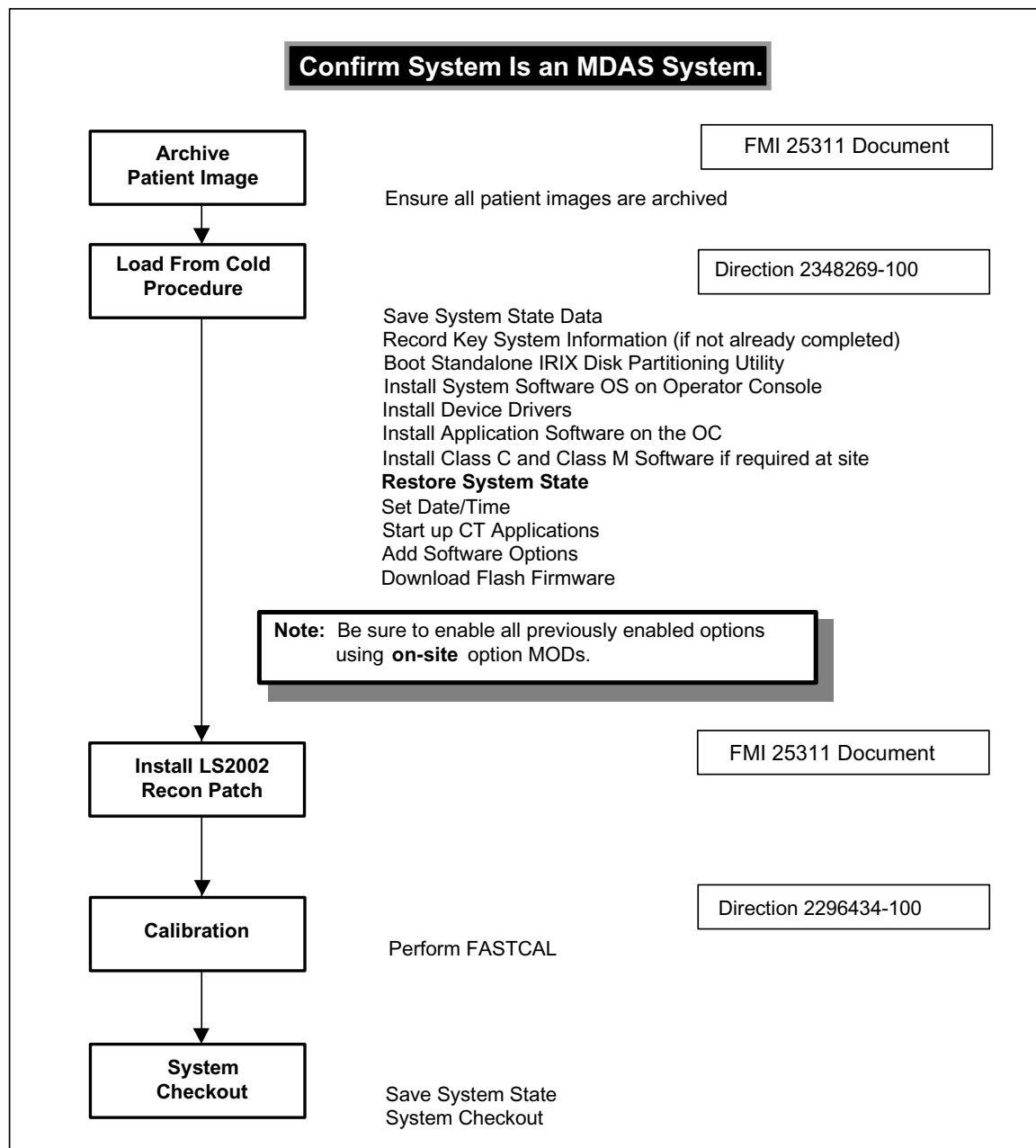

Procedure

1 Archive Patient Images

Verify that the customer has archived all patient data.

2 Software Load Procedures (LFC)

Perform the Load From Cold (LFC) procedure as detailed in Chapter 1 of *LS2002 Software Installation Procedures* (2348269-100).



NOTICE It is required that a “Restore STATE” be performed following the load of Class C and Class M software. To facilitate this, it is recommended that the Class C and Class M be loaded following the Applications Load (Section 3.6).

3 Installing the LS2002 M3 Recon Patch

3.1 Installation

- 1.) Shut down applications by:
Service -> Utilities -> Application Shutdown
- 2.) Open a UNIX shell and become root by:
type: `su -`
- 3.) Enter the root password.
- 4.) Insert the LS2002 M3 Recon Patch CD in the CD-ROM Drive:
 - a.) type: `mountCD /CDROM`
 - b.) type: `/CDROM/install_ls2002_m3_patch`The install messages will appear:
.....
.....
- 5.) When the installation is complete:
 - a.) type: `unmountCD`
 - b.) Reboot the system

3.2 Verification of the Software Patch Installation

To verify Patch is installed in OC:

type: `showprods | grep ls2002_patch_b`

The output will look similar to the following:

```
ls2002_patch_b 07/10/2002 Patch for ls2002 M3
ls2002_patch_b.img 07/10/2002 Patch files for recon
ls2002_patch_b .img.ocroot 07/10/2002 Files that go onto ocroot
```


4 Calibration

Perform FASTCAL by selecting Daily Prep on the bottom of the left monitor, then selecting FASTCAL from the displayed selections.



NOTICE You may get a pop-up indicating: “System detects an alignment or change was performed on the System”. If this happens, press the Ok button and continue with the FASTCAL.

5 System Checkout and Save State

- 1.) Refer to the Software Installation Procedures Manual, and perform a Save System State. This will save the new calibrations as well as the customer and service protocols.
- 2.) Refer to the Software Installation Procedures Manual, and perform the system performance checks as described in the manual. Ensure all networking and filming processes are still functional.
- 3.) Complete and return Red Installation Product Locator Cards.

6 Appendix: Uninstalling the LS2002 M3 Recon Patch

- 1.) Shut down applications by:
Service -> Utilities -> Application Shutdown
- 2.) Open a UNIX shell and become root by:
type: `su -`
Enter the root password.
- 3.) Insert the LS2002 M3 Recon Patch CD in the CDROM Drive:
 - a.) type: `mountCD /CDROM`
 - b.) type: `/CDROM/uninstall_ls2002_m3_patch`The uninstall messages will appear:
.....
.....
- 4.) When the installation is complete:
 - a.) type: `cd`
 - b.) type: `unmountCD`
 - c.) Reboot the system

This Page Intentionally Blank.