

**ADMINISTRATIVE INFORMATION****FMI No: 25319****FMI Title: LightSpeed 16 DAS Upgrade**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS-AM:	Issue Date:
☒ Mandatory	GEMS-E:	Close Date:
☐ Mandatory on Request	GEMS-A:	

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	FMI # To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2369572-100****Do NOT leave unsecured on site.**



Completion Certification Report

Mail To

GEMS - Americas

GE Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe

GEMS - E / ESC
FMI Admin. 322
BP34
F 78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25319 FMI Title: LightSpeed 16 DAS Upgrade

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code:

- 1.) Modification Complete or previously modified.
- 2.) FMI not required per FMI instruction.
- 3.) Not modified. Unit in storage. FMI kit returned to stock.
- 4.) Modification refused or equipment altered. Customer signature and title required.
- 5.) Unit location unknown or unit not in area. Product Locator notified.

**Customer Refusal or Equipment Altered
- Modification Pending.**

Installer
Comments

Effectivity

Affected CT equipment includes only those LightSpeed gantries (P/N 2339985) described below:

Site	FDO	Gantry #	System ID	DAS #	Model #
Stanford	0970-53492	583	415723SCT3	MDAS16-0001	2332633
Charite	0930-48079	1157	CC0618CT02	MDAS16-0007	2332633
FMLH	0970-65139	946	414805CT16	MDAS16-0003	2332633
JFK	0970-82047	937	732632ect	MDAS16-0004	2332633
Shanghai	0930-08530	1091		MDAS16-0005	2332633
Keio	0930-11723	1275		MDAS16-0015	2332633

Purpose

The purpose of this FMI is to replace early versions of the DAS subsystem with current production hardware. FMI 25320 must be performed in parallel with this procedure. There are cable support brackets that need to be installed on the new DAS Chassis before they are mounted on the gantry.

Prerequisite FMIs

FMI Number	FMI Title
FMI 25320	LightSpeed 16 Power Supply Cable Upgrade

Related FMIs

FMI Number	FMI Title
FMI 25321	LightSpeed 16 Power Supply Bracket and Tech Pubs Upgrade
FMI 25322	LightSpeed 16 DAS Duct Upgrade

Furnished Materials

Qty	Description	Part #
1	FMI Document	2369572-100
1	Right DAS Chassis	2332630
1	Center DAS Chassis	2332631
1	Left DAS Chassis	2332632
2	Chassis to chassis cables	2338891
1	20 amp fuse	46-170021p81
1	20 amp label	2371217

Required Tools

Qty	Description	Part #
1	FMI Document	2369572-100
1	Field Torque Wrench Kit	Field Supplied
1	10mm Hex Key, to fit torque wrench	Field Supplied
1	8mm Hex Key, to fit torque wrench	Field Supplied
1	Side cutter	Field Supplied
1	Flat blade screw driver	Field Supplied
1	DAS/Detector Service Kit (FMI 25321)	2344539

People Power

7 hours labor and 1 hour of travel.

Procedure

1 HW Adjustment and Installation

1.1 DAS Chassis (Left, Center and Right)

1.1.1 Chassis Removal (Left, Center or Right)



NOTICE



Ensure that you are properly grounded by using the appropriate ESD wrist strap and cord connected to a good ground point in the Gantry.

- 1.) Move the table cradle to the home position, and lower the table to its lowest elevation.
- 2.) Remove gantry right side cover, and turn OFF all three (3) service switches on the **STC** backplane. Remove the rest of the gantry covers.
- 3.) Turn off all 3 switches (Axial Drive, **HVDC**, 120VAC) on the **STC** backplane.
- 4.) Position the bad chassis at 6 O'clock.
- 5.) Lock the gantry in position, using the rotational lock.
- 6.) Remove the **DAS** air plenum, per FMI 25322 instruction.
- 7.) Put on wrist strap and use good **ESD** precautions.
- 8.) Using Detector Handling precautions (verify use of **ESD** strap, grounding and Nitrile gloves), disconnect Detector flexes from the backplanes using the flex extraction tool found in the 16 Slice DAS/Detector Service kit.
- 9.) Disconnect the power cables:
- 10.) Disconnect the data cables:
- 11.) Remove Center Chassis, four (4) large 12mm cap screws, and remove chassis.
- 12.) Place chassis on an **ESD** pad (to protect converter cards).

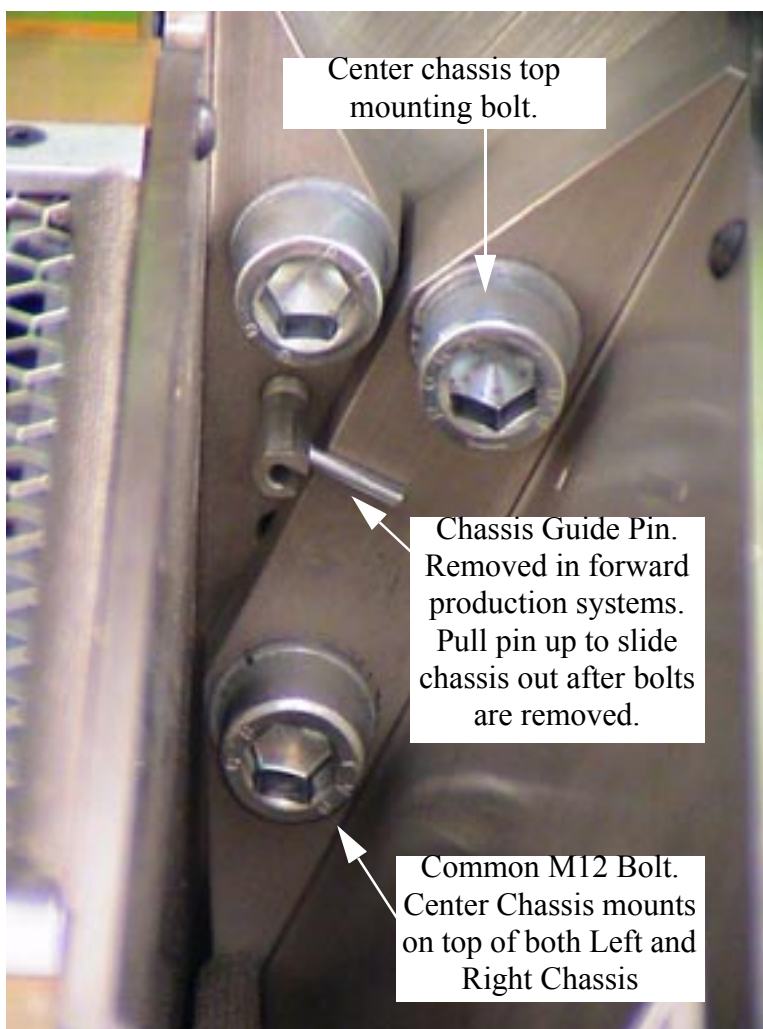


Figure 1 - DAS Chassis mounting bolts, between right and center chassis

- The center chassis shares one bolt each with the left and right chassis.
 - Be careful not to damage detector flexes.
- 13.) Remove right and left chassis using same process as center chassis. Place chassis on an **ESD** pad (to protect converter cards).

1.1.2 Chassis Installation (Left, Center and Right)



NOTICE



Ensure that you are properly grounded by using the appropriate ESD wrist strap and cord connected to a good ground point in the Gantry.

- 1.) Install new DAS chassis. (Center chassis last). Secure loosely with mounting bolts. Use special care:
 - Do not smash or otherwise damage Detector flexes
 - Ensure that ALL Detector flexes are in front of backplane.
- 2.) Secure by using four (4) large 12mm Cap screws, with flat and lock washers. Torque each chassis mounting screw to 66.4 N-m (48.9 lbs-ft).
- 3.) Connect the appropriate power cables:
- 4.) Connect the new chassis to chassis Data Cables.

- Fold each cable, one for the left/ctr and for the right/ctr chassis assemblies. See [Figure 2](#) and [Figure 3](#) for folding details.

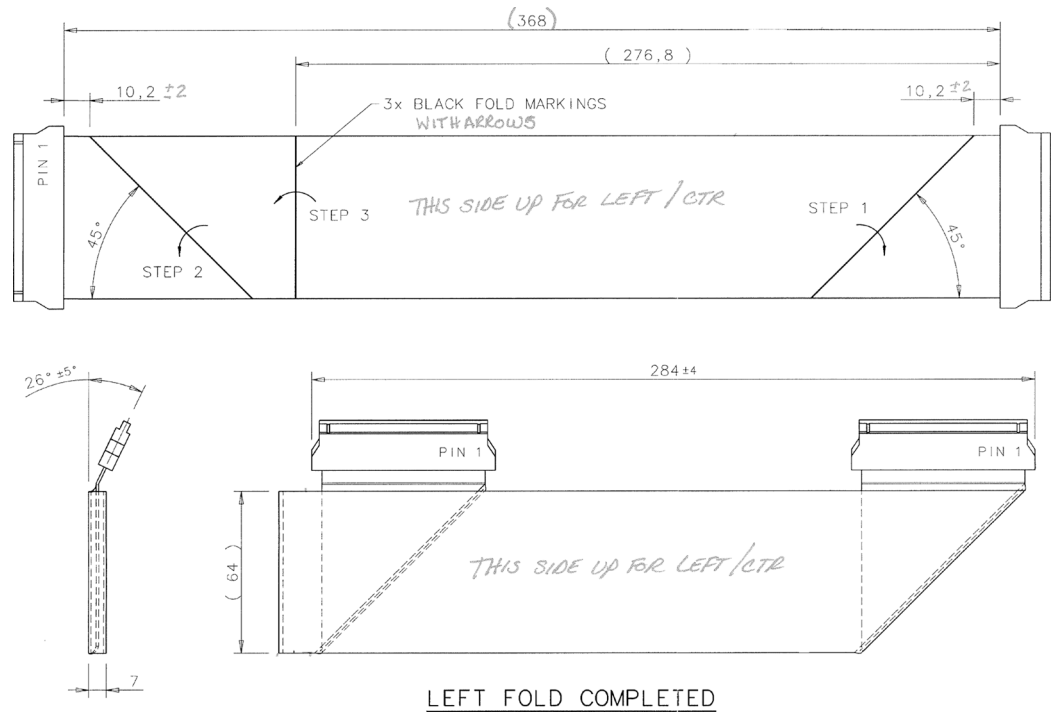


Figure 2 - Left Chassis Fold

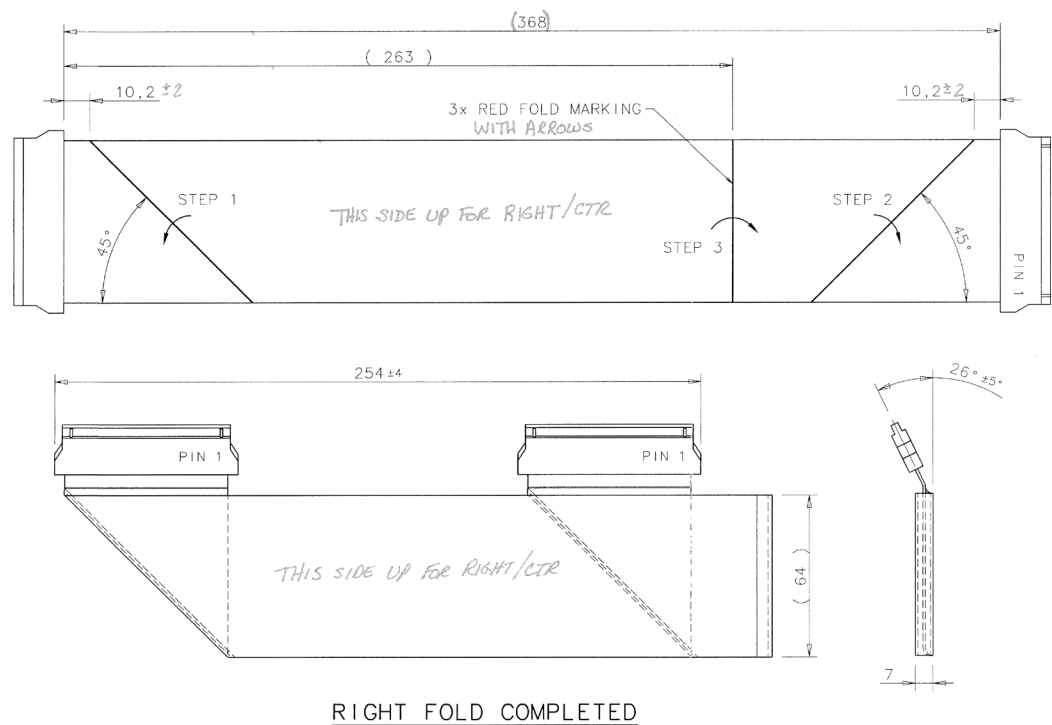


Figure 3 - Right Chassis Fold

- Secure cables in panduit cable clamps.

- 5.) Transfer the Converter boards from the replaced chassis to the new chassis.
 - a.) Confirm **ESD** practices
 - b.) Remove Chassis board cover
 - c.) Remove 1 board at a time and transfer each board from the original chassis to the new chassis in the same board slot location.
- 6.) Install Detector flexes with your Nitrile glove protected fingers.
 - a.) Confirm **ESD** and Detector flex handling practices. Use **ESD** Wrist strap and Nitrile gloves.
 - b.) Interposers and Dust Shields are shipped with each new chassis.
- 7.) Install chassis board cover.

For the right chassis, reattach the fiber-optic cable to the transmit (TX) port.
- 8.) Install plenum, per FMI 25322 instruction.
- 9.) Disengage rotational lock.
- 10.) Replace the "DAS Power" fuse with the provided 20 amp fuse.
- 11.) Affix new fuse label to the fuse box, over the "DAS POWER" silk-screen.
- 12.) Restore power.
- 13.) See Retest Matrix for tests after replacing a **DAS** Chassis.
- 14.) Install / Close all covers.

1.2 Retest Matrix

- 1.) **DAS** Detector Integration
 - a.) X-Ray Verification
 - b.) Microphonics
 - c.) Pop/Noise
- 2.) Air Cal
- 3.) Image Quality Checks
- 4.) Save State

2 Finishing Up

- 1.) Package all old parts and send to the appropriate GE Recycling Center.
- 2.) Complete FMI paperwork, and submit as required.

This Page Left Blank Intentionally.