

**ADMINISTRATIVE INFORMATION****FMI No: 25321****FMI Title: LightSpeed 16 Power Supply Bracket and Tech Pubs Upgrade**

| <b>FMI TYPE</b>                                                                                                                                                  | <b>CHARGE CODES</b>                                         | <b>IMPLEMENTATION DATES</b>                  |
|------------------------------------------------------------------------------------------------------------------------------------------------------------------|-------------------------------------------------------------|----------------------------------------------|
| <input checked="" type="checkbox"/> <b>Mandatory Safety</b><br><input type="checkbox"/> <b>Mandatory</b><br><input type="checkbox"/> <b>Mandatory on Request</b> | <b>GEMS-AM:</b><br><br><b>GEMS-E:</b><br><br><b>GEMS-A:</b> | <b>Issue Date:</b><br><br><b>Close Date:</b> |

|                                                       |                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                          |
|-------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <input checked="" type="checkbox"/> <b>Staged FMI</b> | <b>Customer List</b> provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.<br><br><input checked="" type="checkbox"/> <b>Parts</b> and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.<br><br><input type="checkbox"/> <b>Parts Ordering:</b> GEMS - Americas, contact your Region logistics coordinator<br>GEMS - Europe, Order FMI kits from GPO/EPL Buc<br>GEMS - Asia, Contact your region logistics coordinator<br><br><input type="checkbox"/> <b>No Parts</b> are required to complete this FMI. |
| <input type="checkbox"/> <b>Batch with</b>            | <b>FMI #</b><br>To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                         |
| <input type="checkbox"/> <b>Non-Staged</b>            | <b>No Customer List is provided.</b> Each service location must identify affected units and place part orders as indicated in the Special Instructions.                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                                  |

**Technical Support**

| <b>GEMS - Americas</b>                                                                                                               | <b>GEMS - Europe</b>                                                                                                                                | <b>GEMS - Asia</b>                                                                       |
|--------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------|
| <b>GE Medical Systems</b><br><b>Insite Center</b><br><b>(414) 524-5300</b><br><b>800-321-7973</b><br><b>Milwaukee, WI 53201, USA</b> | <b>GEMS - E / ESC</b><br><b>BP34</b><br><b>F 78533 BUC CEDEX, FRANCE</b><br><b>Tel: (33 1) 39 20 00 07 or ESC</b><br><b>Fax: (33 1) 30 70 99 70</b> | <b>Asia Support Center</b><br><b>Phone: 81-426-56-0033</b><br><b>Fax: 81-426-48-2905</b> |

**Special Instructions:****Unit Tracked:****FMI Instructions Part No.: 2369574-100****Do NOT leave unsecured on site.**



# Completion Certification Report

## Mail To

**GEMS - Americas**

GE Medical Systems  
Attn: FMI Admin. W-523  
P.O. Box 414  
Milwaukee, WI 53201

**GEMS - Europe**

GEMS - E / ESC  
FMI Admin. 322  
BP34  
F 78533 Buc Cedex, France

Send to your  
Country FMI or Service  
Administrator

FMI No.: 25321 FMI Title: LightSpeed 16 Power Supply Bracket and Tech Pubs Upgrade

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

### Unit Tracked:

| Model Number | Serial Number | Compl. Date<br>DD/MM/YY | FMI Comp.<br>Code | Service Hours |        |
|--------------|---------------|-------------------------|-------------------|---------------|--------|
|              |               |                         |                   | Labor         | Travel |
|              |               |                         |                   |               |        |
|              |               |                         |                   |               |        |
|              |               |                         |                   |               |        |
|              |               |                         |                   |               |        |

### FMI Completion Code:

- 1.) Modification Complete or previously modified.
- 2.) FMI not required per FMI instruction.
- 3.) Not modified. Unit in storage. FMI kit returned to stock.
- 4.) Modification refused or equipment altered. Customer signature and title required.
- 5.) Unit location unknown or unit not in area. Product Locator notified.

### Customer Refusal or Equipment Altered - Modification Pending.

Installer  
Comments

## ***Effectivity***

Affected CT Equipment includes the LightSpeed Gantry P/N 2339985 with the following Gantry #s:

| Site                     | FDO                        | Gantry # | System ID  | DAS #       | Model # |
|--------------------------|----------------------------|----------|------------|-------------|---------|
| Stanford                 | 0970-53492                 | 583      | 415723SCT3 | MDAS16-0001 | 2332633 |
| Charite                  | 0930-48079                 | 1157     | CC0618CT02 | MDAS16-0007 | 2332633 |
| FMLH                     | 0970-65139                 | 946      | 414805CT16 | MDAS16-0003 | 2332633 |
| JFK                      | 0970-82047                 | 937      | 732632ect  | MDAS16-0004 | 2332633 |
| Shanghai                 | 0930-08530                 | 1091     |            | MDAS16-0005 | 2332633 |
| U of M                   | 0970-59963                 | 1259     | 734936ct1  | MDAS16-0012 | 2332633 |
| Keio                     | 0930-11723                 | 1275     |            | MDAS16-0015 | 2332633 |
| CCN                      | 0930-23220                 | 1306     | M2322212   | MDAS16-0018 | 2332633 |
| Mayo                     | 0970-13536                 | 1314     | ct81441    | MDAS16-0019 | 2332633 |
| Cardiff                  | 0930-13963                 | 1328     | 00243CTS05 | MDAS16-0020 | 2332633 |
| Kanawanza                | 0930-10356                 | 1265     |            | MDAS16-0014 | 2332633 |
| Bry Sur Marne            | 0930-9932 /<br>0930-43660? | 1333     |            | MDAS16-0025 | 2332633 |
| St Catherine             | 0930-79709                 | 1335     | M2842530   | MDAS16-0024 | 2332633 |
| Cahors                   | 8937-21720                 | 1297     | M4168119   | MDAS16-0017 | 2332633 |
| Rambouillet              | 0930-16974                 | 1348     | M4200607   | MDAS16-0027 | 2332633 |
| Thibodaux                | 0970-64756                 | 1361     | 504447tlsp | MDAS16-0029 | 2332633 |
| L'Orangerie              | 08937-31977                | 1337     | M1455820   | MDAS16-0023 | 2332633 |
| Solingen                 | 0930-046354                | 1381     |            | MDAS16-0030 | 2332633 |
| Groenemeyer              | 0930-056209                | 1382     | UC1590CT01 | MDAS16-0032 | 2332633 |
| Chatswood<br>Diag Centre | 0930-49007                 | 1383     |            | MDAS16-0034 | 2332633 |
| Seirei<br>Hamamatsu      | 0930-032504                | 1384     | QZ0009     | MDAS16-0035 | 2332633 |
| Tokyo Police             | 0930-10360                 | 1386     |            | MDAS16-0037 | 2332633 |
| Nuffield                 | 930-22208                  | 1373     | 00557CTS01 | MDAS16-0036 | 2332633 |
| Skovde                   | 0930-056242                | 1387     | R002CT02   | MDAS16-0040 | 2332633 |
| Medicsight               | 0930-044457                | 1392     |            | MDAS16-0044 | 2332633 |

## ***Purpose***

The purpose of this FMI is to replace early versions of the DAS +/- 5 vdc analog power supply brackets. These early brackets had machining defects which could result in the Service Engineer damaging the detector mounting plate by stripping the M12 female threads. Additionally current service documentation and a DAS/Detector Service Kit are to be given to the customer for support of their LightSpeed 16 CT Scanner.

## ***Furnished Materials***

| Qty | Description                                     | Part #        |
|-----|-------------------------------------------------|---------------|
| 1   | FMI Document                                    | 2369574-100   |
| 4   | Power Supply Bracket                            | 2332623       |
| 6   | screw - hexagon socket 12mm X 90 mm             | 46-312577P105 |
| 6   | washer - lock spring 12.2 mm X 18.1 mm          | 46-313575P7   |
| 6   | washer plain - normal 12.5 mm X 20 mm           | 2275713-8     |
| 16  | screw hex socket 4 mm X 12 mm                   | 46-328417P4   |
| 16  | washer - plain - normal 4.3 mm X 9 mm           | 46-328430P2   |
| 16  | Nut hexagon 4 mm                                | 46-328425P1   |
| 4   | screw countersunk head 4mm 6mm c-sunk flat head | 46-328422P12  |
| 4   | cable tie mount adhesive back                   | 46-208747P2   |
| 1   | Loctite 271 10 cc bottle                        | 46-170684P2   |
| 1   | Documentation Collector                         | 2354016       |
| 1   | 16 Slice DAS/Detector Service Kit               | 2344539       |

## ***Required Tools***

| Qty | Description                        | Part #         |
|-----|------------------------------------|----------------|
| 1   | FMI Document                       | 2369574-100    |
| 1   | Field Torque Wrench Kit            | Field Supplied |
| 1   | 10mm Hex Key, to fit torque wrench | Field Supplied |
| 1   | 3mm Hex Key, to fit torque wrench  | Field Supplied |
| 1   | 7mm Open End Wrench                | Field Supplied |

## ***People Power***

4 hours labor and 1 hour of travel.

## ***Procedure***

### **1 Before You Begin**

Do the following first.

- 1.) Check DAS power supply values using DASTOOLS to establish reference.  
See Direction 2341960-100, Rev 01, Chapter 2, Section 4.5.5.1, for details.
- 2.) Inventory FMI kit to verify all parts are present.

## 2 HW Adjustment and Installation



### DANGER



**POTENTIAL FOR ELECTROCUTION OR SHOCK. HAZARDOUS VOLTAGES ARE PRESENT IN THE GANTRY, EVEN WITH ALL STC SERVICE SWITCHES “OFF”. SERVICING OF THE MDAS 16 SHOULD BE PERFORMED BY QUALIFIED PERSONNEL ONLY.**

Note: It may be easier to access the power supplies if you remove the gantry top and rear covers. Then tilt the gantry all the way back to 30 degrees.

- 1.) Position the cradle in the home position and raise the table to its highest elevation.
- 2.) Remove gantry right side cover, and turn OFF the Axial and HVDC service switches on the **STC** backplane.
- 3.) Remove the gantry left, top and rear covers.

Note: Tilt the gantry 30 degrees forward or backward for easiest access.

- 4.) Perform System Lockout/Tagout procedure.  
See Direction 2341960-100, Rev 01, Chapter 1, Section 4.1, for details.
- 5.) Rotate gantry so that the **DAS** assembly is in a position where the power supply can be easily accessed. This could be anywhere between the 9 o'clock and 3 o'clock positions, depending on which power supply is to be replaced and the availability of a stepping stool (to work from the top of the gantry).



### WARNING



**POTENTIAL FOR PERSONAL INJURY. USE THE ROTATIONAL LOCK TO PREVENT UNCOMMANDED MOVEMENT OF THE GANTRY.**

- 6.) Lock the gantry in position, using the rotational lock




### WARNING



**POWER SUPPLY MAY BE HOT TO TOUCH. USE CAUTION TO PREVENT POSSIBLE BURN INJURY FROM HOT SURFACE CONTACT.**

- 7.) Disconnect the rotating power harness as required.
- 8.) Remove the plastic cover by removing four (4) each: 4mm cap screws, flat and lock washers.
- 9.) Disconnect the AC and DC connectors.
- 10.) Remove the three (3) 12mm cap screws that secure the power supply to the detector plate.  
[Two screws are located at one end of the power supply assembly, and the third is at the other end.] Lift the power supply assembly out of the gantry.
- 11.) Remove the eight (8) 4mm cap screws and associated hardware securing the brackets to the power supply chassis. Discard this hardware.
- 12.) Install the new brackets using the new 4mm x 12mm cap screws, washers and nuts. Apply Loctite 271 to each 4mm x 12mm cap screw. Torque to 2.3 N-m (20.4 lb.-in). [See Figure 1.](#)
- 13.) Install the cable tie mounts using the new 4mm x 6mm countersunk screws. Torque to 1.4 N-m (12.4 lb.-in). [See Figure 1.](#)

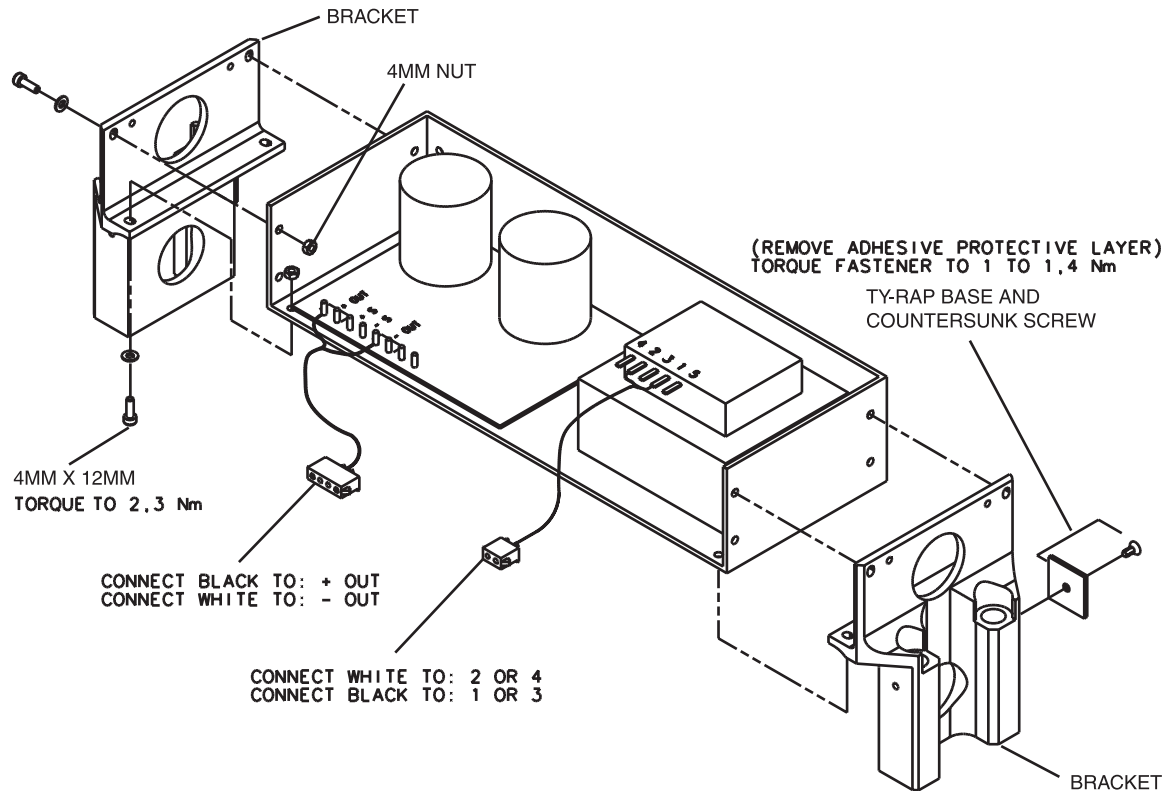

**Figure 1 - 5 VDC Power Supply Mounting Bracket Assembly Drawing**

- 14.) Install the power supply assembly to the detector plate, using the NEW three (3) 12 mm cap screws, one (1) flat and one (1) lock washers. Torque the screws to 66 N-m (48.68 lbs-ft).



**CAUTION**

**Use only one (1) lock washer and one (1) spring washer per 12mm x 90mm cap screw. Use these and only these washers as they set the correct thread engagement for reliable fastener clamping force. Failure to do this can result in loose power supply (12mm cap screw bottoms out) or stripped detector plate female threads (not enough thread engagement).**

- 15.) Reconnect the supply output and supply power.
- 16.) Attach the plastic cover to the assembly, using 4mm cap screws, flat and lock washers. Torque the screws to 0.7 N-m (6.2 lbs-in).
- 17.) Secure the rotating harness cables as before.
- 18.) Restore power to the gantry and verify power supply output voltage. Adjust supply output, if necessary.

### 3 Finishing Up

- 1.) Package all old parts and send to the appropriate GE Recycling Center.  
Deliver to the Customer the updated service documentation and 16 Slice DAS/Detector Service kit for the support of their LightSpeed 16 scanner.
- 2.) Complete FMI paperwork and submit as required.