

**ADMINISTRATIVE INFORMATION****FMI No: 25322****FMI Title: LightSpeed 16 DAS Duct Upgrade**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS-AM:	Issue Date:
☒ Mandatory	GEMS-E:	Close Date:
☐ Mandatory on Request	GEMS-A:	

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	FMI # To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2369575-100****Do NOT leave unsecured on site.**



Completion Certification Report

Mail To

GEMS - Americas

GE Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe

GEMS - E / ESC
FMI Admin. 322
BP34
F 78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25322 FMI Title: LightSpeed 16 DAS Duct Upgrade

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code:

- 1.) Modification Complete or previously modified.
- 2.) FMI not required per FMI instruction.
- 3.) Not modified. Unit in storage. FMI kit returned to stock.
- 4.) Modification refused or equipment altered. Customer signature and title required.
- 5.) Unit location unknown or unit not in area. Product Locator notified.

**Customer Refusal or Equipment Altered
- Modification Pending.**

Installer
Comments

Effectivity

Affected CT Equipment includes the LightSpeed Gantry P/N 2339985 with the following serial #s:

Site	FDO	Gantry #	System ID	DAS #	Model #
Stanford	0970-53492	583	415723SCT3	MDAS16-0001	2332633
Charite	0930-48079	1157	CC0618CT02	MDAS16-0007	2332633
FMLH	0970-65139	946	414805CT16	MDAS16-0003	2332633
JFK	0970-82047	937	732632ect	MDAS16-0004	2332633
Shanghai	0930-08530	1091		MDAS16-0005	2332633
U of M	0970-59963	1259	734936ct1	MDAS16-0012	2332633
Keio	0930-11723	1275		MDAS16-0015	2332633
CCN	0930-23220	1306	M2322212	MDAS16-0018	2332633
Cardiff	0930-13963	1328	00243CTS05	MDAS16-0020	2332633
Kanawanza	0930-10356	1265		MDAS16-0014	2332633
Bry Sur Marne	0930-9932 / 0930-43660?	1333		MDAS16-0025	2332633

Purpose

The purpose of this FMI is to replace the initial version of LightSpeed 16 systems, shipped with a pre-production DAS air plenum, with current production hardware. Please note that the air filters from the old air plenum will not work with the new air plenum.

Furnished Materials

Qty	Description	Part #
1	FMI Document	2369575-100
1	DAS Duct Assembly	2333651

Required Tools

Qty	Description	Part #
1	FMI Document	2369575-100
1	Field Torque Wrench Kit	Field Supplied
1	8mm Hex Key, to fit torque wrench	Field Supplied
1	Side cutter	Field Supplied

People Power

1.5 hours labor and 1 hour of travel.

Procedure

1 HW Adjustment and Installation

1.1 Plenum Removal

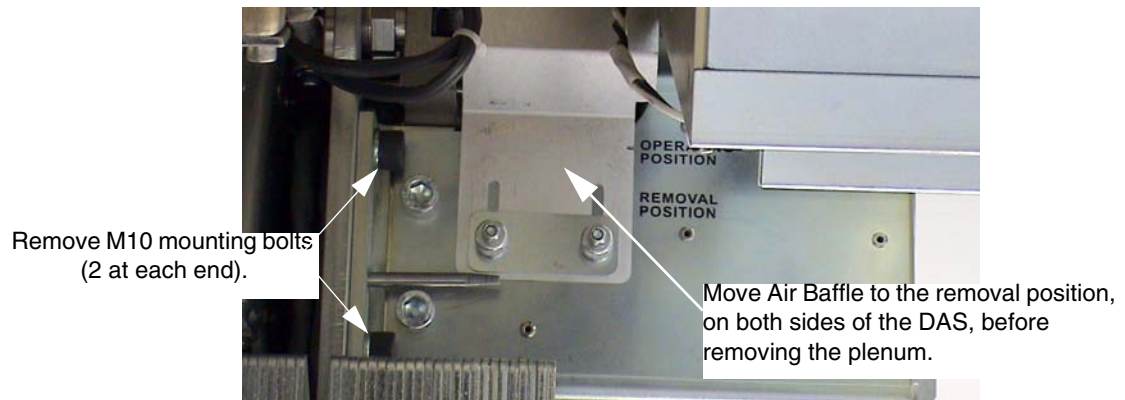


Figure 1 - DAS Air Plenum: Baffle & Mounting Bolts

- 1.) Move the table cradle to the home position, and lower the table to its lowest elevation.
- 2.) Remove gantry side cover, and turn OFF all three (3) service switches on the **STC** backplane.
- 3.) Remove the gantry top, side and front covers.
- 4.) Position the DAS at 12 O'clock.
- 5.) Lock the gantry in position, using the rotational lock.
- 6.) Disconnect Fan AC power at Connector J22.
- 7.) Set air baffles on either end of the duct to the "Removal Position" (see [Figure 1](#)):
 - a.) Loosen knurl nuts
 - b.) Slide baffle into position
 - c.) Hand tighten knurl nuts
- 8.) Remove four (4) M10 mounting bolts (2 at each end). See [Figure 1](#).
- 9.) Carefully slide the plenum out, without catching wires or Detector flexes.
- 10.) Set the plenum aside, on a clean surface.

1.2 Plenum Installation

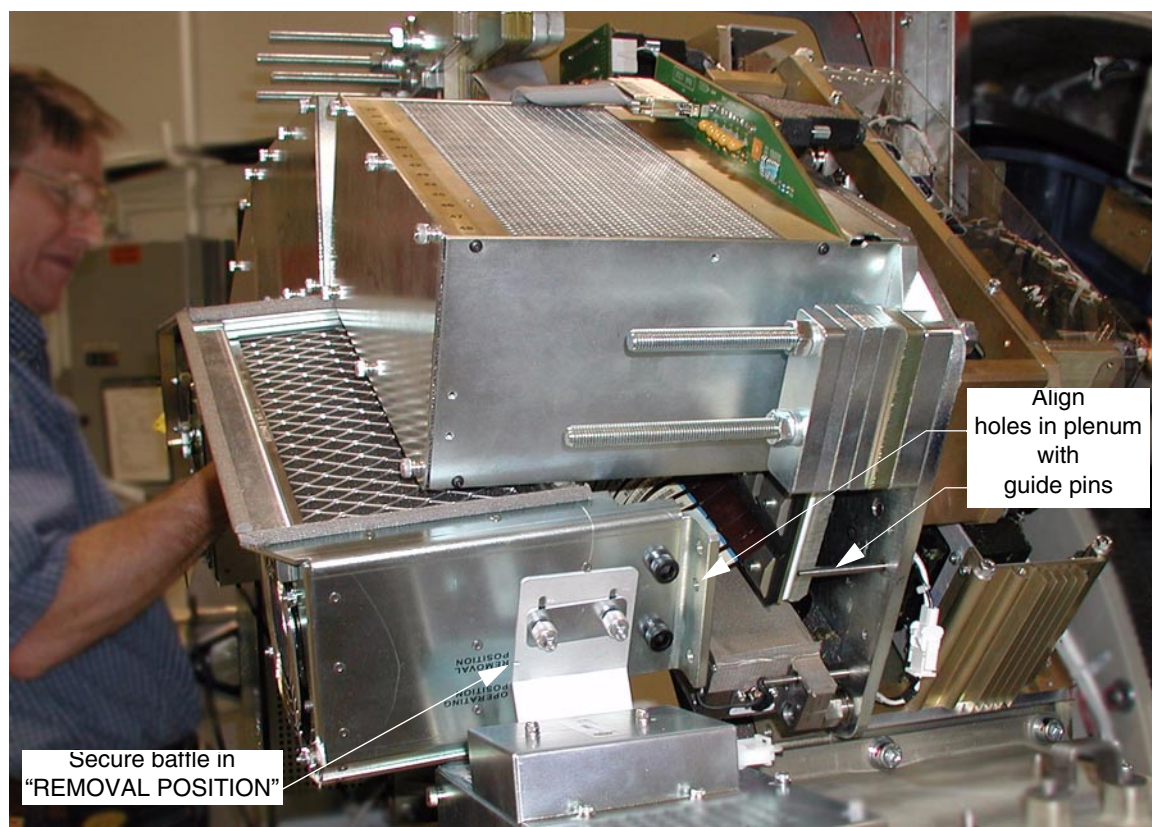


Figure 2 - DAS Air Plenum Installation

- 1.) Ensure that baffles are secured in "REMOVAL POSITION".
- 2.) Lift plenum into place, and align guide holes with guide pins.
- 3.) Slide plenum onto guide pins.
- 4.) Fasten plenum securely into place, using the four (4) M10 mounting bolts and four lock washers. Torque to 38.4 N-m (28.3 lbs-in).
- 5.) Return baffles to the "OPERATING POSITION".
- 6.) Connect Fan AC power at Connector J22.

2 Finishing Up

- 1.) Package all old parts and send to the appropriate GE Recycling Center.
- 2.) Complete FMI paperwork and submit, as required.

This Page Left Blank Intentionally.