

**ADMINISTRATIVE INFORMATION****FMI No: 25332****FMI Title: LS2002 M5 Software Upgrade - LightSpeed 2X SDAS**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support		
GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:

Unit Tracked: 2260619 & 2260619-3

FMI Instructions Part No.: 2399103-100

Do NOT leave unsecured on site

FMI Completion Certification Report

Mail To

GEMS - Americas

G.E. Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe

GEMS - Europe
FMI Admin. 322
BP34
F78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25332

FMI Title: LS2002 M5 Software Upgrade - LightSpeed 2X SDAS

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments

Effectivity

The following CT LightSpeed Equipment is affected:

LightSpeed 2X systems with SDAS. Gantry Model Type Helios 2:

- 2260619 &
- 2260619-3.



NOTICE: Do Not Load on PET/CT Systems, Linux Consoles (GOC 2) or GRE Consoles (GOC 3).

Purpose

The purpose of this FMI is to upgrade LightSpeed IRIX-based CT scanners to the latest version software. This software includes safety and quality fixes.

Safety (PSR / FPR) fixes include:

- PSR 1286436 – Auto mA table values,
- PSR 1286629 – Auto mA noise index,
- PSR 13000489 – Scan abort with prep delay, and
- PSR 1286829 – Wrong Patient Information.

Quality fixes include:

- Scout not displayed,
- ring artifact,
- remote tilt,
- align light,
- Autolink page slow,
- Batch Install Image not displayed,
- Image out of order,
- Slow image-biopsy,
- Transfer to Pacs Error,
- Service Desktop Error,
- autoMA-concerns (table),
- Modem,
- autoMA-noise (GUI),
- death of stclf,

- axial not ready,
- Smart Prep fails to start,
- .rhosts file corrupted,
- scan data missing in retro,
- Square VP not saved,
- AutoArchive,
- peripheral gate hang,
- ali pacs (mekesson),
- Can't recover MOD images,
- not accept image - Phillips,
- Show Local lines delay, and
- Head IQ.

New delivered features are:

- AutomA GUI,
- Repeat Series,
- Re-Start / Shutdown Procedure, and
- Anonymous Scan Data Save.

Furnished Materials

Qty	Part number	Description
1	2399103-100	FMI document
1	2395178	LS2002 M5 SDAS Software Collector containing:
1	2395182	LS2002 M5 Application CD Collector (3 CDs)
1	2356012	LightSpeed OS Software Collector (3 CDs)
1	46-302200P2	Product Locator
1	2348269-100	LS2002 Software Install Document
1	2243321-100	LightSpeed 2X Apps Tips and Workarounds
1	2369740-200	CT Lightspeed Series Learning and Reference Guide


Tools Required

Qty	Description	Part #
1	LightSpeed Class C Software (Tab Distribution)	2395172
1	LightSpeed Class M Software (Tab Distribution)	2395173
1	Variviewer Software Option (Customer Site)	2347142
1	Auto mA Software Option (Customer Site)	2356009
1	Power 440 Software Option (Customer Site)	2345411
-	Customer Purchased Options (Customer Site)	-

People Power

5 hours labor and 1 hours of travel

FMI Installation Sequence Guide

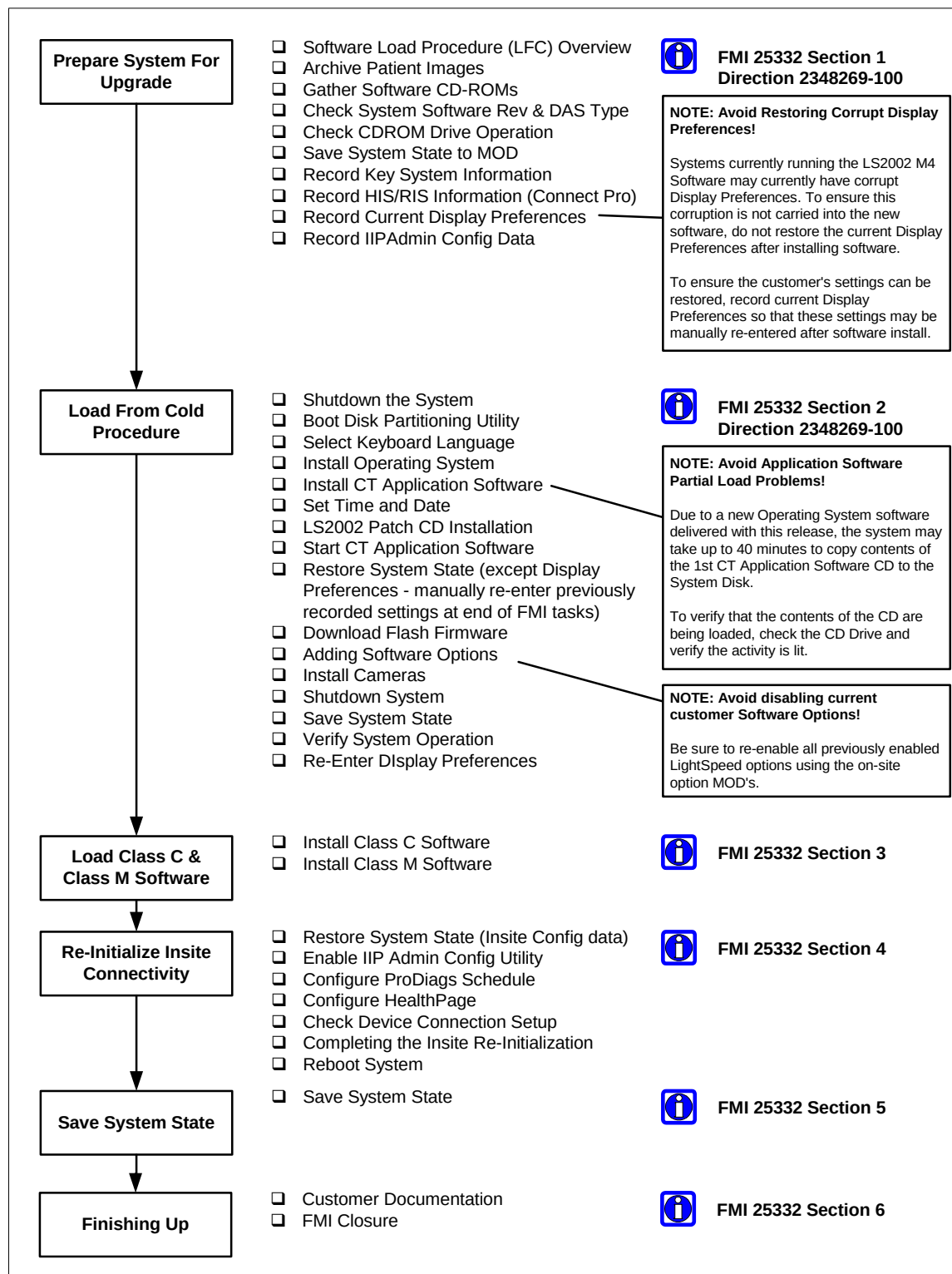

Procedure

Introduction

This FMI includes the following procedures:

- Preparing the System for Software Upgrade
- Load From Cold

Section 1

Prepare System for Software Upgrade

Follow steps in S/W Load Procedure 2348269-100 Chapter 1, Section 2 Before You Begin:

- Check System Software Version (Version 2.1, 2.2 or 2.3) & Verify DAS Type (SDAS)
- Gather Software CD-ROMs
- Check CD-ROM Drive Operation
- Save System State to MOD
- Record KEY System Information
- Record IP address for Broadband (Only if the system has Broadband)

BroadBand Insite Connected

If your system is not broadband connected to Insite, omit this section.

If your system has been Insite Checked Out as a broadband system, you will need to re-enter the Insite Gateway IP Address 2399103-100ss into the IIP Admin Configuration utility when reconfiguring Insite after the Load From Cold process.

To ensure that you have entered the correct IP Address of the Insite dedicated Gateway, confirm this IP Address with one of the following resources:

- ask the Hospital Network Administrator or IT contact.
- in the United States, check for the IP Address at the following web site or by calling the regular checkout number (800-321-7937):
<http://gein2.med.ge.com/vpni/secure/index.jsp>
- in Europe, contact the regular toll free number to reach the OLCE.

Perform this investigation prior to arriving on site.

- Record IP address below:

- IP Address ____ . ____ . ____ . ____

Archive Image Data

The following procedures (Software Load From Cold) will destroy all image data on the system. Prior to starting the software load, make sure that the customer has reconstructed and archived all images before proceeding.

System Operation Verification

SETUP AND PERFORM AT LEAST ONE AXIAL SCAN TO VERIFY THAT THE SYSTEM OPERATED PROPERLY BEFORE BEGINNING THE UPGRADE.

Insite Connectivity Info (found in /usr/g/insite .insiteINFO)

If your system is not connected to Insite, omit this section.

Check CT LightSpeed Insite Connectivity Settings as follows:

1. Open a shell
2. Type:
cd /usr/g/insite <ENTER>
more .insiteINFO <ENTER>

From the usr/g/insite .insiteINFO file, record the following Insite Connectivity configuration parameters:

Note: If you have Broadband be sure to record the BroadBand Insite Gateway IP for recording later. If you don't have BroadBand you can skip this step

Parameter	Line Name in /usr/g/insite .insiteINFO	Record Parameter Value	Notes
Connection Type	INS_ConnectionType		
Modem Type	INS_ModemType		Modem Connectivity Only
Modem Phone Number	INS_PhoneNumber		Modem Connectivity Only
Modem Dial Prefix	INS_DialPrefix		Modem Connectivity Only
Modem Tone or Pulse Setting	INS_ToneOrPulse		Modem Connectivity Only
Modem Country	INS_Country		Modem Connectivity Only
Modem String	INS_ModemString		Modem Connectivity Only
Modem TTY Port	INS_TtyPort		Modem Connectivity Only
Modem TTY Speed	INS_TtySpeed		Modem Connectivity Only
BroadBand Insite Gateway IP	INS_GatewayIP	-----	BroadBand Connectivity Only
BroadBand Gateway Mask	INS_GatewayMASK	-----	BroadBand Connectivity Only

INSTALLED OPTIONS (FOUND USING SWOKINSTALL -P)

Check CT LightSpeed Options Settings as follows:

- Open a shell.
- Type: **swokinstall -p <ENTER>**

From the swokinstall -p output, record which of the following LightSpeed system options are installed on the system:

Line Name in output file of swokinstall -p	Enabled (Y/N)	Permanent	Flex Trial	Standard or Purchased Option
Smart Prep				Standard
Smart Speed				Configuration Dependent
Helical Tilt				Standard
Direct-3D				Configuration Dependent
Recon Plus				Configuration Dependent
6.3 MHU				Standard
Thin Twin				Standard
VariViewer				Standard
AutomA				Standard
Power440				Standard
Connect Pro				Purchased Option
DentaScan				Purchased Option
SmartScore Pro				Purchased Option
CardIQ SnapShot				Purchased Option
SmartStep				Purchased Option
CTPerfusion2				Purchased Option
VolumeViewer				Purchased Option
AdvVesselAnalysis				Purchased Option
CTColono				Purchased Option
CardIQ2				Purchased Option

NOTE: Make sure that for each Option enabled, there is an on-site MOD that contains the key(s) for the Options presently enabled. Without the MOD keys, you will not be able to restore options on the system that were previously enabled.

Section 2

Load From Cold Procedure

Software Load Procedure

After system operation has been verified, all images have been archived, System State saved, and a record of key system configuration parameters have been completed, you are ready to upgrade the system software.

Using the LS2002 Software Installation Procedure document (**2348269-100**) and the software, both included in the FMI kit, perform a software load.

Reference: Direction **2348269-100 Chapter 1, Section 3.0 Load Procedure.**

Re-enter Display Preferences.

Section 3

Load Class C and Class M Software

Load Class C Software

Load provided Class C software as follows:

1. Open a shell
2. Enter root as superuser:

At the ctuser prompt, type: **su – <ENTER>**

At the Password prompt, type: **#bigguy <ENTER>**

3. Insert the Class C software CD Drive. Wait for the disk read to STOP.
4. Install the Class C software on the system.

At the root prompt, type: **install.classC <ENTER>**

5. Verify Class C is properly loaded on the system.

At the root prompt, type: **showprods | grep CLASSC**

The shell should display the Class C version as follows:
1104I.2_H4.0M5CLASSC

6. Remove the Class C software CD from Drive.

Load Class M Software

Load provided Class M software as follows:

1. Open a shell
2. Enter root as superuser:

At the ctuser prompt, type: **su – <ENTER>**

At the Password prompt, type: **#bigguy <ENTER>**

3. Insert the Class M software CD Drive. Wait for the disk read to STOP.
4. Install the Class M software on the system.

At the root prompt, type: **install.classM <ENTER>**

5. Verify Class M is properly loaded on the system.

At the root prompt, type: **showprods | grep CLASSM**

The shell should display the Class M version as follows:

1204I.2_H4.0M5CLASSM

6. Remove the Class M software CD from Drive.

Section 4

Re-Initializing Insite Connectivity

This release of software contains a new IIP release (IIP 3.2). Changes in this IIP Release primarily in the following areas:

- Updates to support Multi-Tech modem firmware revision 8.19M
- Removal of IIP_DialPrefix ProDiags task
- Replacement of current iipHealthPage ProDiags task with a new Insite Data Mining task
- Auto Route to AutoSC in Level 5 setup

In addition, the Network Configuration Tab in the Insite Interactive Administration utility has been modified to eliminate the Use Default Gateway selection.

To ensure these IIP 3.2 features are correctly set, an Insite Checkout is NOT required. Because the basic Insite configuration information on the system should be correct from the original Insite Checkout, in most cases, it is only necessary to run ConfigLnk via exiting the IIP Admin Configuration utility Insite Checkout Tab Screen.

Restore System State To Retrieve Class M Configuration Parameters

Once Class M software has been loaded, you must repeat a Restore System State to recover all previous Insite and Class M configuration parameters.

Reference: Direction **2348269-100**, LS2002 Software Installation Procedure, Chapter 1, Section 3.9 Restore System State

NOTE: Make sure that the Key System Information recorded in Section 1 (Insite Connectivity Info) is correctly set as recorded. If not, use the following to:

- IIP Admin Configuration utility to reset Network Configuration Info as described in the following procedures.

Enable Insite Interactive Platform Configuration Utility

See Figure 2, Insite Interactive Platform Configuration Utility Gaining Access and Selecting ProDiags Utility.

Enable the Insite Interactive Platform Configuration utility as follows:

1. Open a shell

2. Enter root as superuser:

At the ctuser prompt, type: **su – <ENTER>**

At the Password prompt, type: **#bigguy <ENTER>**

3. Launch the IIP Admin Configuration utility:

At the root prompt, type: **iipadmin config <ENTER>**

The system displays the Insite Interactive Platform Configuration utility with the Proprietary Agreement Acceptance Screen.

4. Obtain access to the IIP Admin Configuration utility:

On the IIP Admin Configuration utility Proprietary Agreement Acceptance Screen: Click on **Accept** three times.

The IIP Admin Configuration utility displays the Information Tab.

5. Launch the independent ProDiags utility from the IIP Admin Configuration Utility:

On the IIP Admin Configuration utility Information Screen:
Click on the **ProDiags** Tab.

The IIP Admin Configuration utility displays the ProDiags Tab.

On the IIP Admin Configuration utility ProDiags Screen:
Click on the **Custom** button.

The IIP Admin Configuration utility launches the independent ProDiags utility as shown in Figure 3.


Figure 2
Insite Interactive Platform Configuration Utility
Gaining Access and Selecting ProDiags Utility



Proprietary Agreement Screen

1.

Click **Accept** Three Times to gain access to the Insite Interactive Platform Configuration Utility

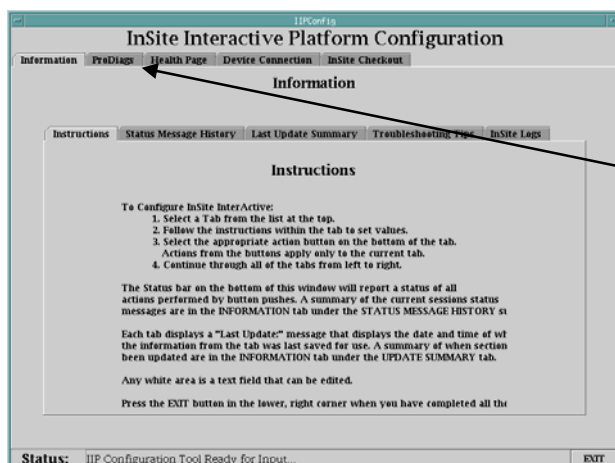

Information Screen

2.

Click on the **ProDiags** Tab to display the IIP Admin Config Utility ProDiags Screen



ProDiags Screen

3.

Click on the **Custom** button to launch the independent ProDiags Utility

Configure ProDiags Schedule

With this release of software, the IIP 3.2 package contains a new data collection package that pushes gesyslog data and a basic system configuration file to the On Line Center RDA Database. This IDM package contains two new ProDiags tasks that must be added to the ProDiags schedule.

The Restore State process (after loading Class M software) restored the older ProDiags schedule which had a similar data collection package called iipHealthPage. This older task must be removed from the schedule.

In addition, there are default set ProDiags tasks in the ProDiags schedule that need to be removed based on On Line Center analysis of AutoSC traffic and usability of information provided by these tasks.

See Figure 3, Insite Interactive Platform Configuration Utility Gaining Access and Selecting ProDiags Utility.

Deleting ProDiags Scheduled Tasks

On the ProDiags utility Schedule Tab, delete the following tasks from the schedule as shown in Figure 3:

FLOG

ChkLogs (all six instances)

IipHealthPage

Use the following sequence for each task deletion:

- Click on the task.
- Click on the **Delete** button at the bottom of the Schedule Tab.
- Click on the **OK** button in the “Are you sure you want to delete” pop-up.

When completed with all deletions, click on the **Apply** button at the bottom of the Schedule Tab. On the “Apply change confirmation” pop-up, click on **Yes**.

Figure 3
ProDiags Graphical Utility
Schedule Tab – Deleting Tasks

ProDiags Utility Schedule Screen

Deleting Tasks From the ProDiags Schedule

1. Click on task to highlight
2. Click on **Delete** button
3. Click on **OK** on the Delete pop-up.
4. Repeat steps 1 through 3 for each task to delete.
5. Click on **Apply** button to apply changes to ProDiags schedule
6. Click on **OK** on the Apply pop-up

Schedule Management

Schedule Management provides the following functionalities

1. Add a task to the schedule
2. Delete a task from the schedule
3. View the schedule

Task	Iterations	Day	Hour	Minute	Type	Time Slot
FLOG	1	Daily	05 AM	00	Background	
PassChg	1	Daily	05 AM	30	Background	
pd_House...	1	Daily	04 AM	30	Background	
pd_ASC_N...	1	Daily	Hourly	28	Background	
ChkLogs	1	Daily	12 PM	00	Background	
ChkLogs	1	Daily	04 PM	00	Background	
ChkLogs	1	Daily	08 PM	00	Background	
ChkLogs	1	Daily	12 AM	00	Background	
ChkLogs	1	Daily	04 AM	00	Background	
ChkLogs	1	Daily	08 AM	00	Background	
healthog	1	Daily	04 AM	45	Background	
IIP_HouseK...	1	Daily	06 AM	12	Background	
IIP_HouseK...	1	Daily	12 PM	12	Background	
IIP_HouseK...	1	Daily	06 PM	12	Background	
IIP_HouseK...	1	Daily	12 AM	12	Background	
idm_send...	1	Daily	07 AM	00	Background	
idm...	1	Daily	Hourly	00	Background	

New Delete Reload Apply Help

Status : Reading data... Done

Delete the following ProDiags tasks from the Schedule:

- FLOG
- ChkLogs (all six instances)
- iipHealthPage

Adding ProDiags Scheduled Tasks

On the ProDiags utility Schedule Tab, add the following tasks to the schedule as shown in Figure 4:

Idm

Idm_send_data

Use the following sequence for the **idm** and **idm_send_data** task additions:

- Click on the **New** button at the bottom of the Schedule Tab.
- From the “Adding New Schedule” pop-up, select the **idm** task and click on **OK**.
- On the “Add Confirmation” pop-up, click on **Yes**.
- Click **OK** on the pop-up requesting an edit of the newly created scheduled item.

Leave the default setting of the **idm** task (Daily, Hourly).

- Click on the **New** button at the bottom of the Schedule Tab.
- From the “Adding New Schedule” pop-up, select the **idm_send_data** task and click on **OK**.
- On the “Add Confirmation” pop-up, click on **Yes**.
- Click **OK** on the pop-up requesting an edit of the newly created scheduled item.
- Click on the “Day” module of the **idm_send_data** task and, using the pull down list, select **Daily**.
- Click on the “Hour” module of the **idm_send_data** task and, using the pull down list, select a time where the system is not in prime use but is not regularly turned off.

When completed with the **idm** and **idm_send_data** additions, click on the **Apply** button at the bottom of the the Schedule Tab. On the “Apply change confirmation” pop-up, click on **Yes**.

Close the ProDiags utility browser by selecting **Close** from the upper left browser window pull down menu. The ProDiags browser should close and the Insite Interactive Platform Configuration utility should display the HealthPage Tab as shown in Figure 5.

Figure 4
ProDiags Graphical Utility
Schedule Tab – Adding New Schedule Tasks

Close the ProDiags utility by selecting **Close** in the pull down menu

ProDiags Utility Schedule Screen

Adding Tasks to the ProDiags Schedule

1. Click on **New** button. Utility displays "Adding New Schedule" pop-up.
2. On the "Adding New Schedule" pop-up, click on the **idm** task to highlight and click on **OK**.
3. On the Add Confirmation pop-up, click **Yes**.
4. On the pop-up requesting a schedule edit of the task, click **OK**.
Leave the default setting of the **idm** task - Daily and Hourly.
5. Repeat Steps 1 through 3 for the **idm_send_data** task.
6. On the ProDiags Schedule Tab, update the **idm_send_data** task schedule as follows:

Day: **Daily**
Hour: **07 AM**
Minute: **00**
7. Click on **Apply** button to apply changes to ProDiags schedule
8. Click on **OK** on the Apply pop-up

Schedule Management

Schedule Management provides the following functionalities

1. Add a task to the schedule
2. Delete a task from the schedule
3. View the schedule

Task	Iterations	Day	Hour	Minute	Type	Time Slot
FLOG	1	Daily	05 AM	00	Background	
PassChg	1	Daily	05 AM	30	Background	
pd_House...	1	Daily	04 AM	30	Background	
pd_ASC_N...	1	Daily	Hourly	28	Background	
ChkLogs	1	Daily	12 PM	00	Background	
ChkLogs	1	Daily	04 PM	00	Background	
ChkLogs	1	Daily	08 PM	00	Background	
ChkLogs	1	Daily	12 AM	00	Background	
ChkLogs	1	Daily	04 AM	00	Background	
ChkLogs	1	Daily	08 AM	00	Background	
healthog	1	Daily	04 AM	45	Background	
IIP_HouseK...	1	Daily	06 AM	12	Background	
IIP_HouseK...	1	Daily	12 PM	12	Background	
IIP_HouseK...	1	Daily	06 PM	12	Background	
IIP_HouseK...	1	Daily	12 AM	12	Background	
idm_send...	1	Daily	07 AM	00	Background	
idm	1	Daily	Hourly	00	Background	

Adding New Schedule

Java Applet Window

Task	Task Description	Type	Scheduled
FLOG	Report Failed Lo...	Background	No
PassChg	Report system p...	Background	Yes
idm	Runs IDM Engine	Background	Yes
idm_send_data	Sends data colle...	Background	Yes
ctrld	Tar and Send th...	Background	No
slprune	Prunes the log r...	Background	No

Add the following ProDiags tasks to the Schedule:

Task:	idm	idm_send_data
Iterations:	1	1
Day:	Daily	Daily
Hour:	Hourly	07 AM
Minute:	00	00
Type:	Background	Background

Configure HealthPage

After completing the ProDiags setup, the IIP Admin Configuration utility automatically switches to the HealthPage Tab as shown in Figure 5.


The Restore State process (after loading Class M software) restored the following information for the HealthPage Reports:

- All e-mail addresses of targeted recipients that were previously loaded
- HealthPage report contents that was previously configured
- HealthPage report schedule that was previously configured

Unless you have recent changes to any of these three items, you can skip the HealthPage Tab of the IIP Admin Configuration utility and select the Device Connection Tab. If you have changes that you want to make (e-mail addresses, report content, report schedule), follow the instructions on the IIP Admin Configuration utility.

Figure 5

Insite Interactive Platform Configuration Utility
HealthPage Tab

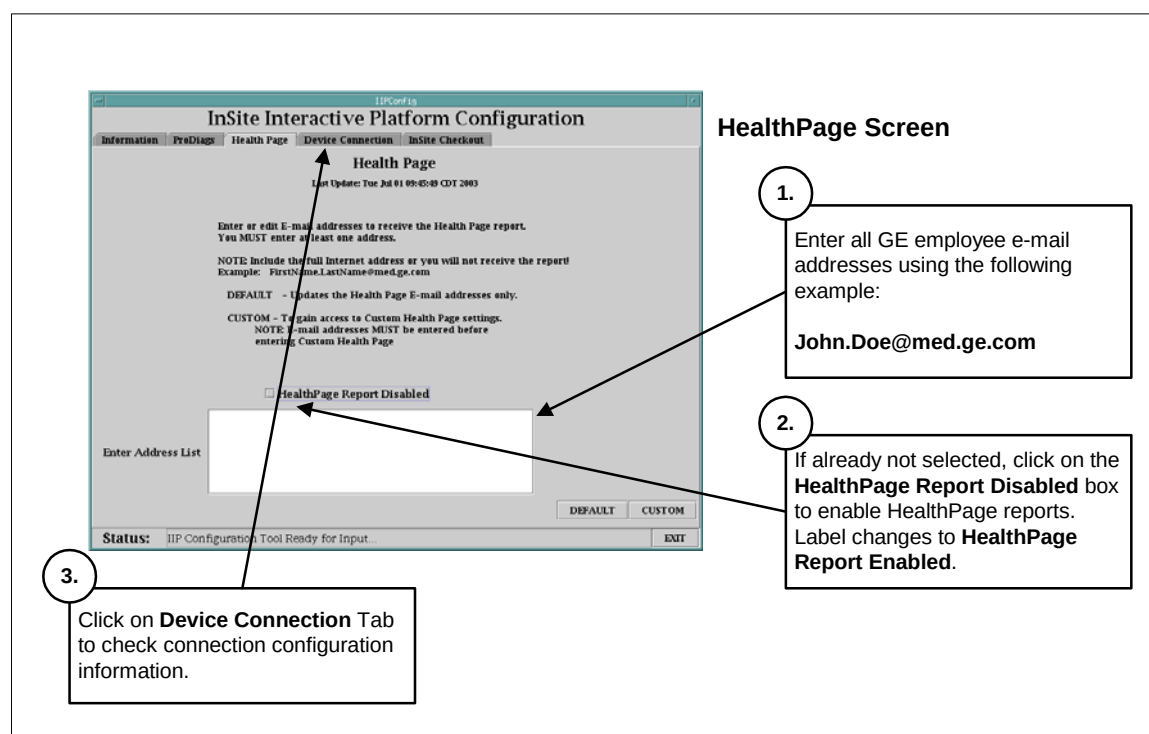

Check Modem Device Connection Setup

When selecting the Device Connection Tab, the IIP Admin Configuration utility automatically displays the Modem Setup Screen as shown in Figure 6.

If your system was configured for Modem connectivity, the Restore State process (after loading Class M software) restored all the Modem setup parameters. To verify, refer to Section 1 of this FMI where you recorded Insite Connectivity Info.

- If there were no Modem parameters in the .insiteINFO file but there were BroadBand parameters, skip to the Check BroadBand Device Connection Setup procedure.
- If there were Modem parameters in the .insiteINFO file but the parameters displayed in the Modem Setup Screen do not match what you previously recorded, enter the correct values as required.

NOTE: The Modem configuration parameters should be pre-loaded from the Restore State process. Do not change settings unless these values were not properly restored.

If required, follow the instructions in Figure 5. The default values are usually the correct values or most commonly used. Parameter selections that do not require customization are:

CPU Serial Port Name
CPU Serial Port Speed

Figure 6**Insite Interactive Platform Configuration Utility
Device Connection Tab – Modem Setup**

Device Connection Modem Screen

1. Select the Dial-out Prefix value recorded from the previous .insiteINFO file (INS_DialPrefix) from the pull down menu. This value is the number you must use to dial out of the hospital facility.

2. Select the Modem Type value recorded from the previous .insiteINFO file (INS_ModemType) from the pull down menu.

3. Select the Dialing Mode value recorded from the previous .insiteINFO file (INS_ToneOrPulse) from the pull down menu.

4. Select the Country value recorded from the previous .insiteINFO file (INS_Country) from the pull down menu.

5. Select the CPU Serial Port Speed value recorded from the previous .insiteINFO file (INS_TtySpeed) from the pull down menu.

6. Select the CPU Serial Port Name value recorded from the previous .insiteINFO file (INS_TtyPort) from the pull down menu.

7. Click on **Apply** to set configuration parameters and display the Insite Checkout Tab.

Connection Device Configuration
Last Update: Tue Jul 01 09:45:58 CDT 2003
Device Connection Type: Modem

Select the InSite Connection Type above, or select the modem settings given below and Press
APPLY - To update the modem and modem configuration files with model and country defaults.
CUSTOM - For other modem settings.

Dial-out Prefix: 9,1
Modem Type: MultiTech MT5634ZBA V.90
CPU Serial Port Name: /dev/ttyS0

Dialing mode: Tone
Country: Default - All Others
CPU Serial Port Speed: 38400

APPLY CUSTOM
Status: Reading Modem Configuration Files...

Check BroadBand Device Connection Setup

If your system was previously BroadBand Insite connected, click on the Network button of the IIP Admin Configuration utility Device Connection Tab to display the Network Setup Screen as shown in Figure 7.

If your system was configured for BroadBand connectivity, the Restore State process (after loading Class M software) restored all the BroadBand setup parameters. To verify, refer to Section 1 of this FMI where you recorded Insite Connectivity Info.


NOTE: The BroadBand Network configuration parameters should be pre-loaded from the Restore State process. However, the on-site FE who performed the original Insite Checkout of the system may have used the hospital's imaging LAN default gateway IP Address (selection was available in the previous release of IIP).

DO NOT USE THE HOSPITAL WORKING LAN DEFAULT GATEWAY IP ADDRESS AS THE INSITE GATEWAY IP ADDRESS OTHER THAN AS A TEMPORARY MEANS TO CONNECT TO THE ON LINE CENTER!!

To find the correct Insite Gateway IP Address, use the contacts listed in Section 1, BroadBand Insite Connected Systems.

Figure 7

Insite Interactive Platform Configuration Utility
Device Connection Tab – Network Setup

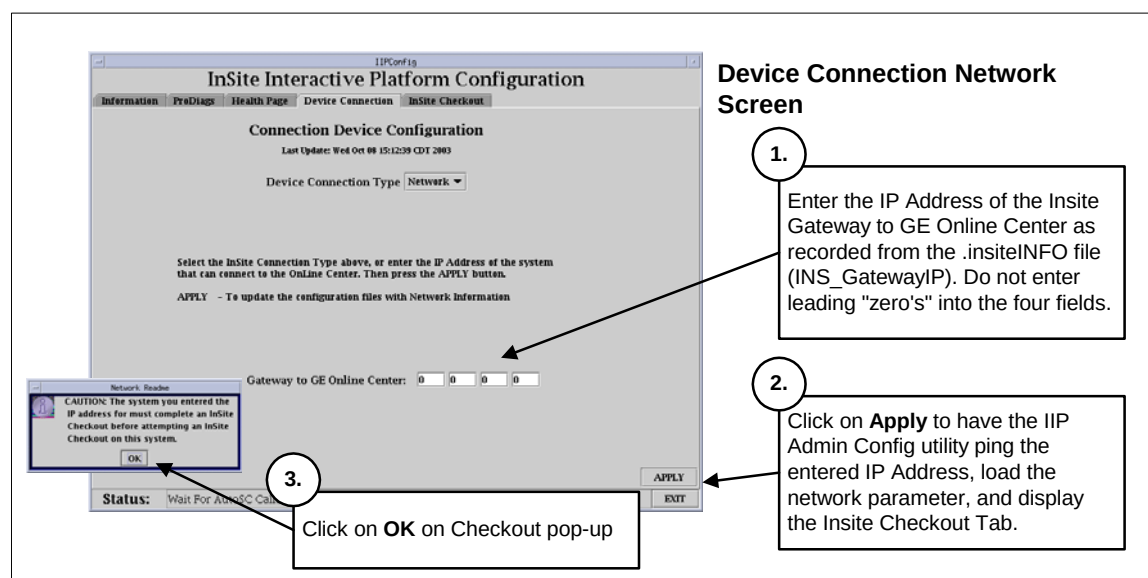

Completing the Insite Re-Initialization


If your system was previously BroadBand Insite connected, click on the Network button of the IIP Admin Configuration utility Device Connection Tab to display the Network Setup Screen as shown in Figure 8.

Follow the instructions provided in Figure 8 to complete the Insite Re-Initialization.

Figure 8

Insite Interactive Platform Configuration Utility

Insite Checkout Tab – Exiting the IIP Admin Configuration Utility

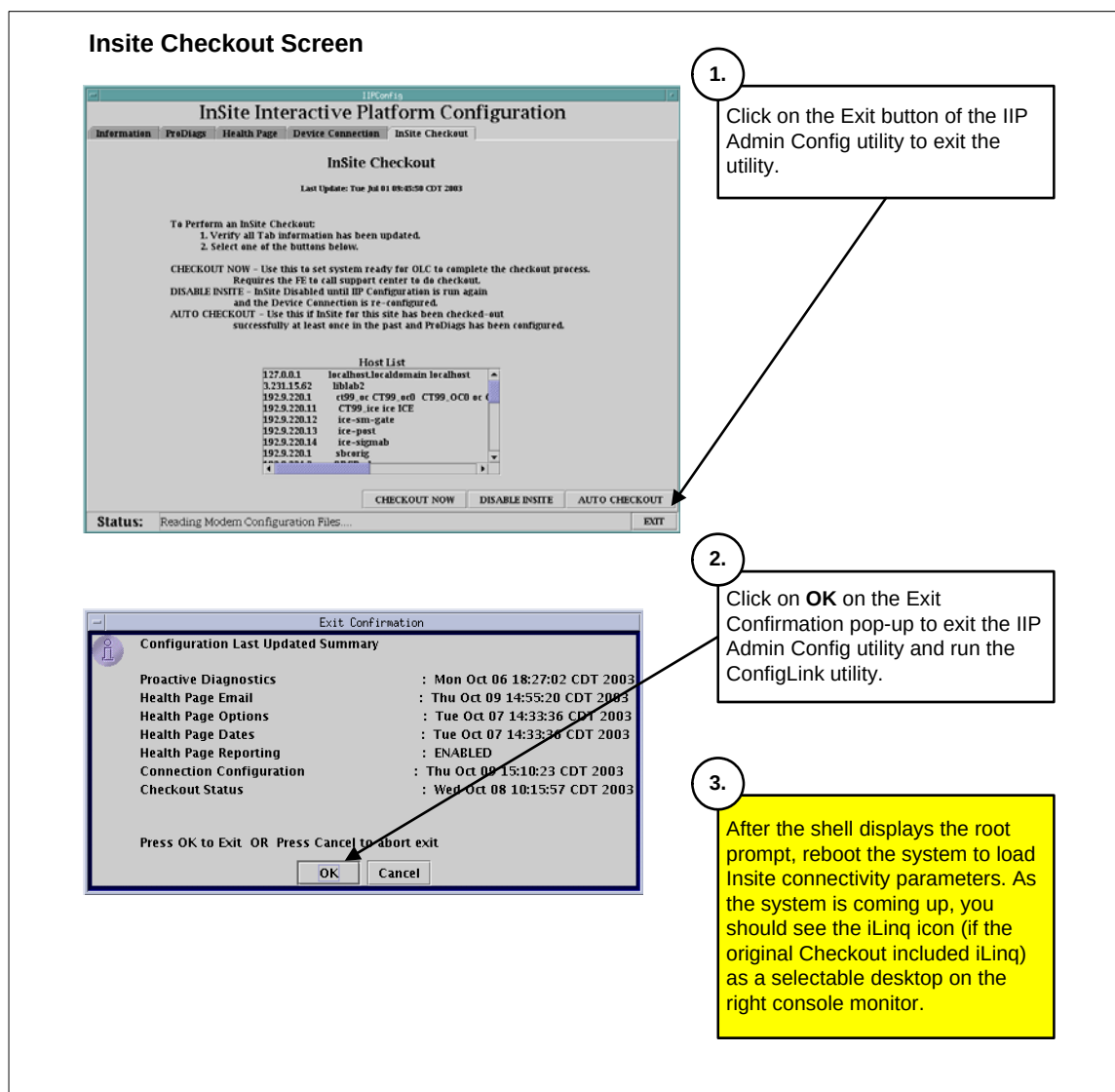

Run the ConfigLink script to populate all configuration information required to fully restore connectivity parameters on the system. Run ConfigLink as follows:

1. 1. Open a shell
2. 2. Enter root as superuser:

At the ctuser prompt, type: `su - <ENTER>`

At the Password prompt, type: `#bigguy <ENTER>`

3. 3. Launch the IIP Admin Configuration utility:

At the root prompt, type: `cd /usr/g/insite <ENTER>`

At the root prompt, type: `ConfigLink <ENTER>`

Wait for the shell to return the root prompt to ensure the ConfigLink process completes. When complete, exit root and close the shell

Section 5

Save System State

At the completion of the Load From Cold Process and Insite Re-Initialization process, ensure all system state parameters are saved using the Save System State Utility.

Reference: **2348269-100** , LS2002 Software Installation Procedure, Chapter 1, Section 2.4 Save System State to MOD.

Finishing Up

1. Customer Documentation

Provide to the **Customer** the Class A Operator/Service documentation and S/W CD-ROM's from this FMI:

Qty	Part number	Description
1	2395178	LS2002 M5 SDAS Software Collector containing:
1	2395182	LS2002 M5 Application CD Collector (3 CDs)
1	2356012	LightSpeed OS Software Collector (3 CDs)
1	2348269-100	LS2002 Software Install Document
1	2243321-100	LightSpeed 2X Apps Tips and Workarounds
1	2369740-200	CT Lightspeed Series Learning and Reference Guide

Review the Applications Tips and Workarounds with the customer.

2. FMI Closure

Close FMI per regional procedure. Note the FMI cannot be closed until Product Locator card complete and GIB updated with all tracked items in this FMI.

Remove FMI document from customer site.