

**ADMINISTRATIVE INFORMATION****FMI No: 25359****FMI Title: LightSpeed / HighSpeed Xstream Console Trackball**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☒ Mandatory on Request	GEMS - A:	Close Date:

☐ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☐	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☒ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support		
GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:** 2377708, -4, -7, -10, -13, -14, -15, 2381796, 5111750**FMI Instructions Part No.: 5130477-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25359

FMI Title: LightSpeed / HighSpeed Xtream Console Trackball

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments

Effectivity

This optional FMI upgrades the trackball on CT LightSpeed, CT HighSpeed, and PET/CT Discovery ST host computer consoles at sites with any one of the following Xstream host computer consoles:

- 2381796 HighSpeed (1 IG) Xstream host computer console assembly
- 5111750 Xstream host computer console assembly for Value 4 Slice, 1 IG, China
- 2377708 LightSpeed 16 Slice (1 IG) Xstream host computer console assembly
- 2377708-4 Xstream host computer console assembly, 1 IG, Hino
- 2377708-7 LightSpeed 4 & 8 Slice (1 IG) Xstream host computer console assembly
- 2377708-10 Xstream host computer console assembly, 4 & 8 Slice, 1 IG, Hino
- 2377708-13 LightSpeed 4 Slice (1 IG) SDAS Xstream host computer console assembly
- 2377708-14 LightSpeed 4 & 8 Slice (1 IG) Xstream host computer console assembly, Israel
- 2377708-15 LightSpeed 16 Slice (1 IG) Xstream host computer console assembly, Israel

Purpose

This optional FMI is to upgrade the Xstream console's PS/2 trackball to an EMC (electromagnetic compatible), USB version trackball to prevent possible console lockup due to ESD (electrostatic discharge).

Furnished Materials

Qty	Description	--FMI Kit 25359, p/n 5130930--	Part #
1	FMI Document		5130477-100
1	Trackball, EMC version		2399753
1	Bracket, cable cover, dual PS/2F + single USB-AF		5123480
1	Cable, USB Trackball		5123478
6	Screw, Panhead, 2.5mm x 5mm		46-328421P3

Tool Required

Screwdriver, Phillips, small blade
Screwdriver, Phillips, medium blade


People Power

Travel: 1 hr

Labor: 1 hr


Procedure

1. Shutdown and switch power OFF at the console.
2. Remove rear console cover.
3. Disconnect all three PS/2 devices from the bracket at the top, rear of the console (trackball, keyboard, and mouse).

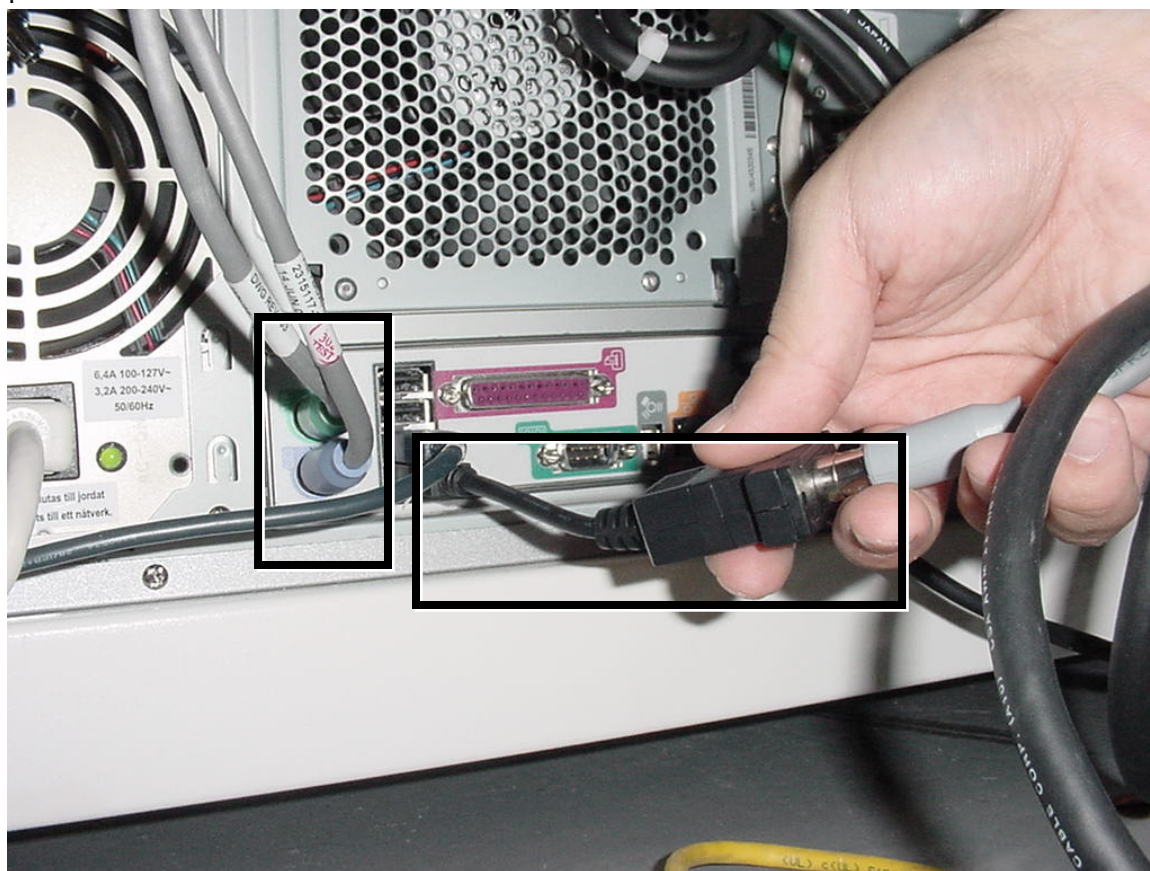

Picture 1: Disconnect Cables from Host PC Computer

4. Disconnect the trackball's PS/2-USB adapter cable, and mouse and keyboard cables from the rear of the console (host PC computer). Refer to the boxed items highlighted in Picture 1.

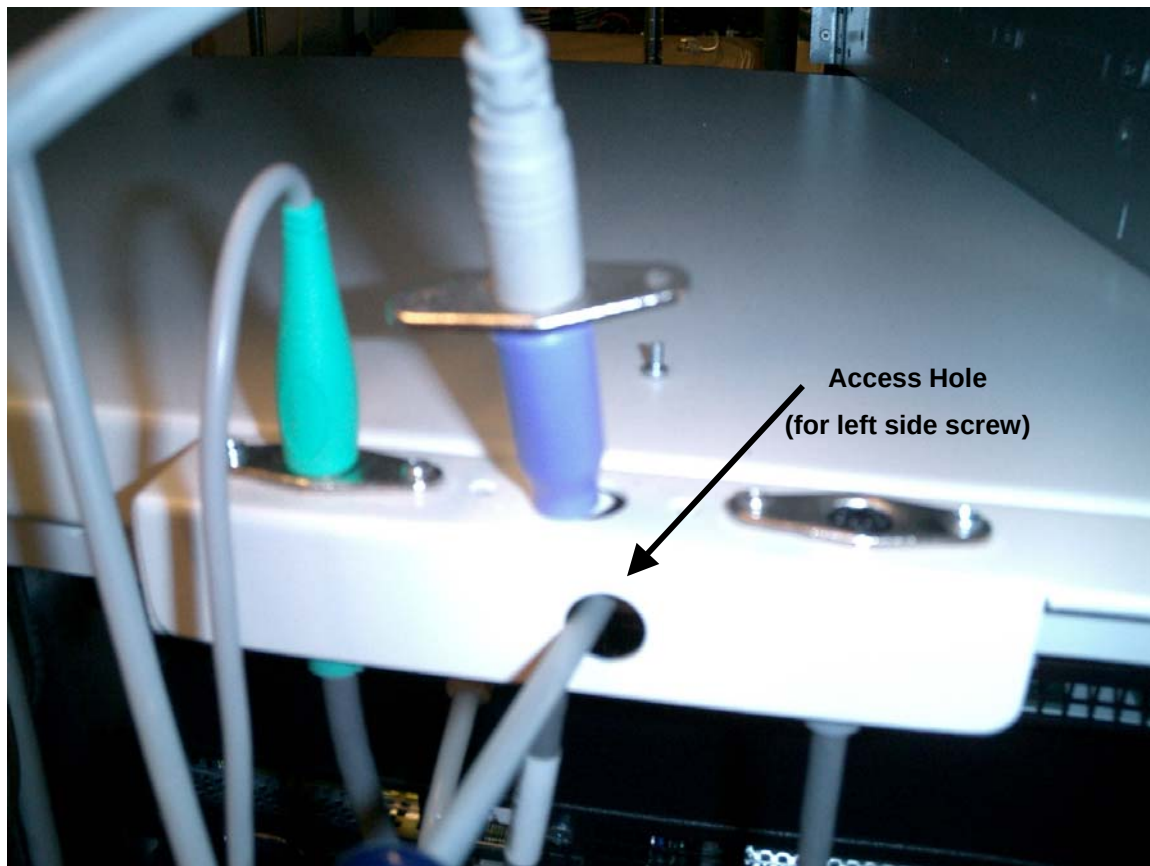


Picture 2: Remove Bracket's Left Side Screw

5. Remove all cable bulkhead connector screws and remove the cables from the bracket by pulling them through their respective holes.
6. Retain the mouse and keyboard cables; they will be reused. Discard the "old" trackball cable; it will be replaced.
7. Loosen bracket's leftmost screw through the access hole; remove it. Refer to Picture 2.



Picture 3: Remove Bracket's Right Side Screw

8. Remove the bracket's rightmost mounting screw. Refer to Picture 3.
9. Attach the new bracket (p/n 5123480, included in kit) to the rear of the console.
Secure it to the console by using the two original screws removed in Step 9 & 10.

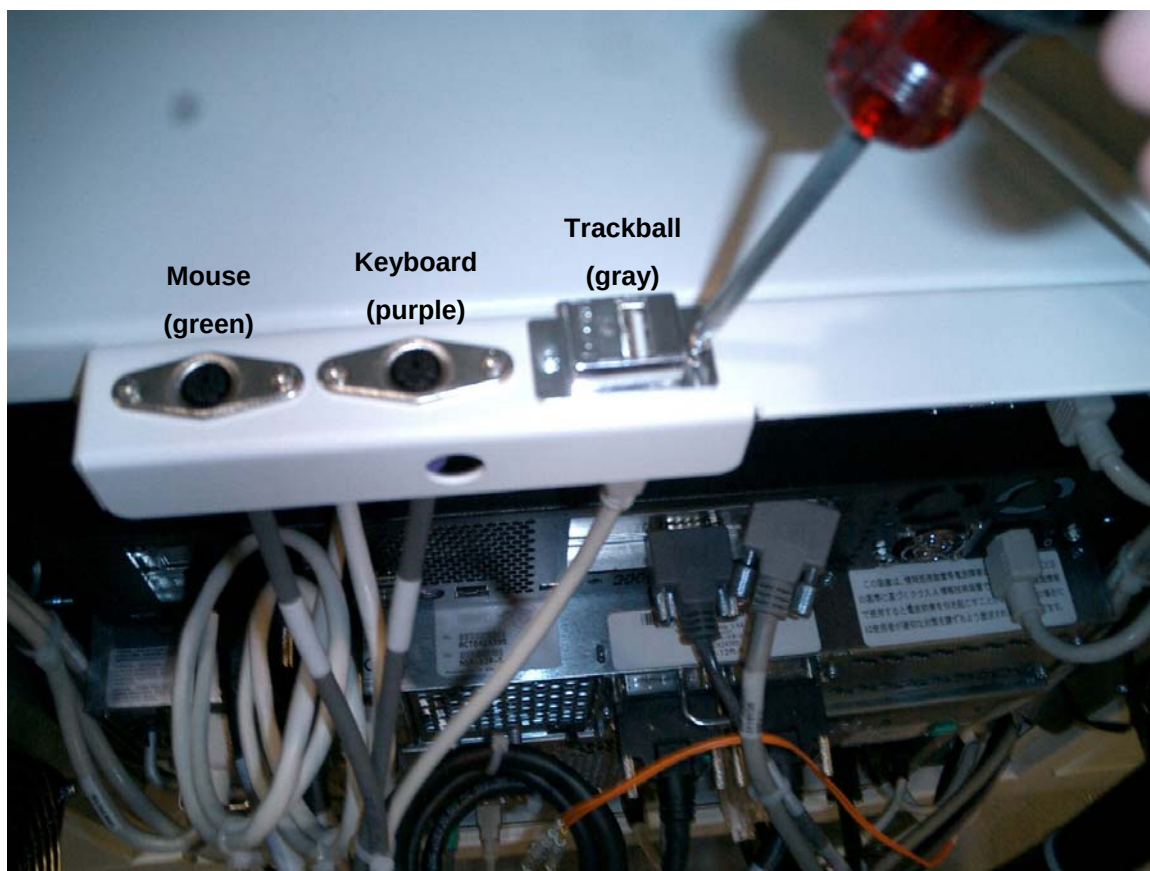


Picture 4: Connect Bulkhead Connectors to Bracket

10. Assemble the two original PS/2 (mouse and keyboard) cables and one new USB (trackball, p/n 5123478, included in kit) cable to the new bracket, p/n 5123480, as illustrated in Picture 4.

Feed each cable through the bracket from the bracket's top side, assembling them to the bracket in the order as illustrated in Picture 4 (green = leftmost, purple=middle, gray=rightmost).

Secure connectors to bracket using six Phillips screws, p/n 46-328421P3, included in kit.

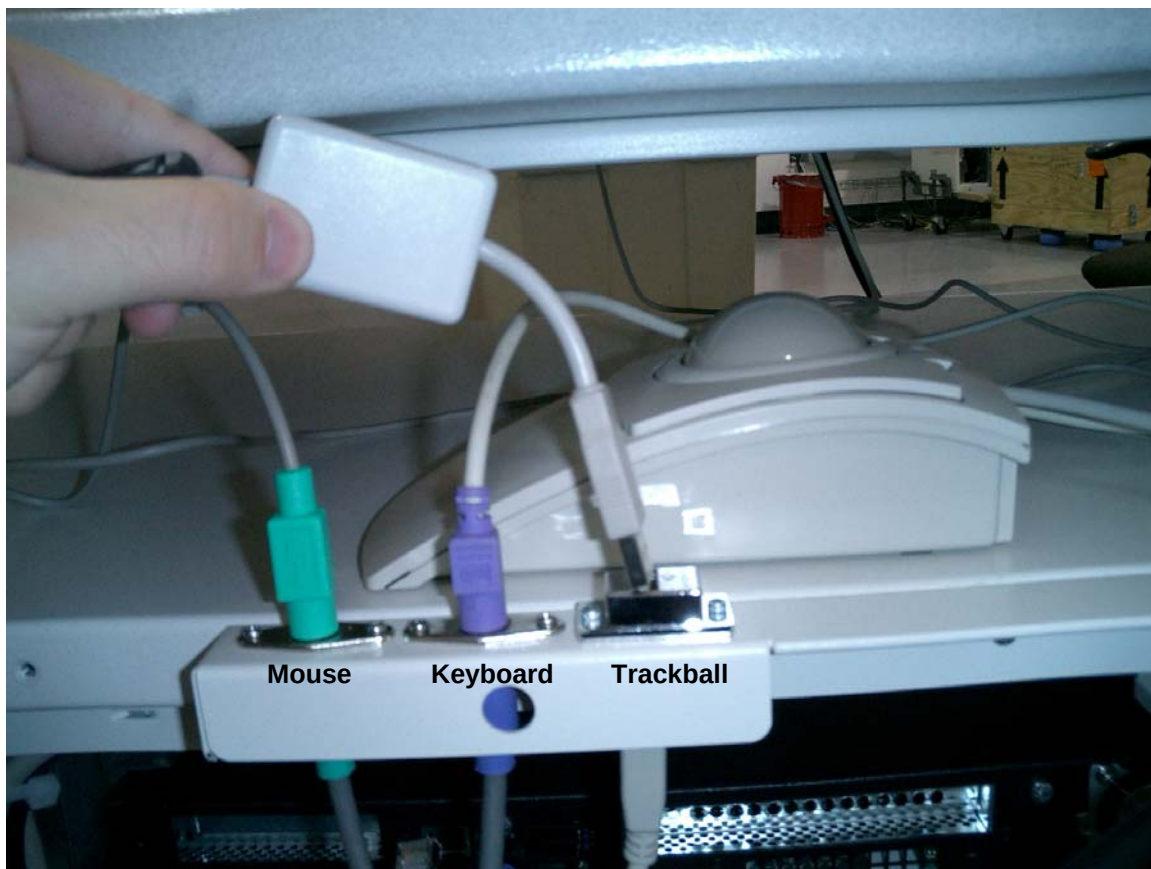

Picture 5: Connect the Mouse, Keyboard, & Trackball

11. Re-connect the mouse, keyboard, and EMC-version trackball (p/n 2399753, included in kit) devices to their respective bulkhead connectors.

Refer to Picture 5.



Picture 6: Connect Cables From Bulkhead to Host PC Computer

12. Connect three cables from the bulkhead / bracket to the host PC computer as illustrated in Picture 6:
 - Mouse (green) to PC's mouse (green) female connector
 - Keyboard (purple) to PC's keyboard (purple) female connector
 - Trackball (gray) to PC's USB port where "old" USB adaptor was connected
13. Re-attach the rear console cover.
14. Turn ON the console's power switch and verify keyboard, mouse, and trackball are functioning.

6 Finishing Up

Do not leave this FMI instruction at the site.

Close FMI per regional procedure.