

**ADMINISTRATIVE INFORMATION****FMI No: 25363****FMI Title: SmarTube Phase 2 SW Patch for Irix/Pegasus systems**

| FMI TYPE                                                    | CHARGE CODES      | IMPLEMENTATION DATES |
|-------------------------------------------------------------|-------------------|----------------------|
| <input checked="" type="checkbox"/> <b>Mandatory Safety</b> | <b>GEMS - AM:</b> | <b>Issue Date:</b>   |
| <input type="checkbox"/> <b>Mandatory</b>                   | <b>GEMS - E:</b>  |                      |
| <input type="checkbox"/> <b>Mandatory on Request</b>        | <b>GEMS - A:</b>  | <b>Close Date:</b>   |

|                                                       |                                                                                                                                                                                                       |
|-------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <input checked="" type="checkbox"/> <b>Staged FMI</b> | <b>Customer List</b> provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole. |
| <input type="checkbox"/>                              | <b>Parts</b> and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.                                                |
| <input type="checkbox"/>                              | <b>Parts Ordering:</b> GEMS - Americas, Contact your Region logistics coordinator<br>GEMS - Europe, Order FMI kits from GPO/EPL Buc<br>GEMS - Asia, Contact your region logistics coordinator         |
| <input type="checkbox"/>                              | <b>No Parts</b> are required to complete this FMI.                                                                                                                                                    |
| <input type="checkbox"/> <b>Batch with</b>            | _____                                                                                                                                                                                                 |
|                                                       | To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).                                     |
| <input type="checkbox"/> <b>Non-Staged</b>            | <b>No Customer List is provided.</b> Each service location must identify affected units, and place part orders as indicated in the Special Instructions.                                              |

**Technical Support**

| GEMS - Americas                                                                                   | GEMS - Europe                                                                                                    | GEMS - Asia                                                         |
|---------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------|
| GE Medical Systems<br>Insite Center<br>(414) 524-5300<br>800-321-7973<br>Milwaukee, WI 53201, USA | GEMS - E / ESC<br>BP34<br>F 78533 BUC CEDEX, FRANCE<br>Tel: (33 1) 39 20 00 07 or ESC<br>Fax: (33 1) 30 70 99 70 | Asia Support Center<br>Phone: 81-426-56-0033<br>Fax: 81-426-48-2905 |

**Special Instructions:****This FMI is for CT Irix/Pegasus (GOC1) systems only. Do not install in CT/PET.****This replaces F.2 software patch. Software prerequisite:**

LS2002M5 S/W LS 2X SDAS systems 2.4.2\_H2.4M5 (FMI 25332, kit # 2399173)

LS 3X MDAS systems 308I.2\_H3.1M5 (FMI 25331, kit # 2399175)

LS 4X MDAS systems 405I.2\_H4.0M5 (FMI 25333, kit # 2399176) Or 405I.3\_H4.0M5.1

HS QXI systems 101I.2HSQXIM5 (FMI 25902, kit # 2401902)

**Unit Tracked:****FMI Instructions Part No.: 5142331-100 Rev 1****Do NOT leave unsecured on site**

# **FMI** Completion Certification Report

## Mail To

**GEMS - Americas**

**GEMS - Europe**

G.E. Medical Systems  
Attn: FMI Admin. W-523  
P.O. Box 414  
Milwaukee, WI 53201

GEMS - Europe  
FMI Admin. 322  
BP34  
F78533 Buc Cedex, France

Send to your  
Country FMI or Service  
Administrator

FMI No.: 25363

FMI Title: SmarTube Phase 2 SW Patch for Irix/Pegasus systems

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

| Model Number | Serial Number | Compl. Date<br>DD/MM/YY | FMI Comp.<br>Code | Service Hours |        |
|--------------|---------------|-------------------------|-------------------|---------------|--------|
|              |               |                         |                   | Labor         | Travel |
|              |               |                         |                   |               |        |
|              |               |                         |                   |               |        |
|              |               |                         |                   |               |        |
|              |               |                         |                   |               |        |

### FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not

### Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments

## ***Effectivity***

CT systems using Irix/Pegasus. Console numbers to be tracked:

- 2266832, 2266832-2, 2266832-3, 2266832-4
- 2304732, 2304732-2, 2304732-3, 2304732-4
- 2341104, 2341104-2, 2341104-3, 2341104-4
- 2341923, 2341923-2

***Note: This FMI is not applicable to CT/PET Systems.***

## ***Purpose***

The purpose of this FMI is to load Patch G.6 on IRIX/Pegasus (GOC1) systems running LS2002 M5 software. This patch includes all the fixes and enhancements of LS2002 M5 patch F.2 (optional FMI 25353).

The software includes SmarTube Phase 2 features (see improvements below).

## ***Improvements***

The F.2 patch fixes shown below are included in the new G.6 patch.


### **Product Safety Problems Fixed:**

- Scout Image Fails to recon  
This was caused by characters added at the end of the Exam Description when the shift or control key was pressed along with the enter key.
- Additional slices added when scans are prescribed with gantry tilt  
Improvements have been made in Graphic Rx and ScanRx
- Show Localizer lines are in wrong location if user skips series quickly
- DLP is higher than expected caused by selecting Repeat Series Quickly.  
The dose reported was incorrect

**Other problems fixed:**

- SmartPrep Diagnostic delay not as expected.

The diagnostic delay is adjusted at Confirm to account for table movement required between the Monitor and Scan phase locations. If the Diagnostic Delay prescribed is less than what the system can do, the expected delay will be posted at Confirm as shown below:

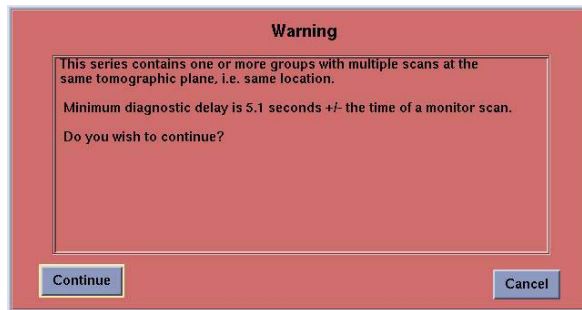

- SmartPrep fails to start
- Improvements have been made in communication between the ScanRx and Display Processes to assure reliable use of SmartPrep
- Start Locations have odd numbers.  
Adjustments to Start Locations from Graphic Rx will be rounded off in .25mm increments. This will give more meaningful scan locations for the user.
- Images can't be retrieved from PACS
- System shuts down with a ScanRx crash
- Dcserver crashes with Storage Commitment requests
- Images are filmed cut off  
This was caused when filming images with GSE from a non-square viewport.
- Update Apache HTTP server to be more secure
- Clean up Pro Diags so Reli Monitor does not fail

**LS2002M5 Software Patch G.6:****Enhancements:****SmarTube**

There are further enhancements to the SmarTube feature in this release. Additional software has been added to make a selection for tube identification based on tube type,

1. GE-Supplied x-ray tube (Performix™).
2. All other tube types.

If a Non-GE-Supplied tube is found, a warning dialog box will be posted at system start-up (see Appendix C).

Advisory messages will be posted on systems with Non-GE-Supplied tubes during Fast Cal and on the Dose Report screen save image indicating that information may not be totally accurate because the full characteristics of the tube may not match those of a GE-Supplied x-ray tube. If the installed x-ray tube is a GE-Supplied tube but the field service engineer did not have all the appropriate information to complete verification, a message (see Appendix C) will be displayed for 30 days or until the field service engineer is able to complete tube installation verification within the 30 day window.

If the verification is not completed in the 30 day window, software will revert to unrecognized tube state and messages for non-validated x-ray tube will be seen during FastCal and on the Dose Report screen save image.

**Dose Report Screen Save: (Octane Dual Processor Systems only)**

Dose Report Screen Save feature has been implemented in this release. The Dose Report is a screen save image that is automatically saved to Series 999 at normal End Exam. The Dose Report is a copy of the information displayed in the Dose Information Area on the View Edit screen. The report includes the CTDIvol, Dose Length Product (DLP) and total DLP for the exam. Scouts will be listed, however, there is no recognized method for measuring dose for scouts so the information for these series will always be blank. If there is more than one group in a series, the series will be listed multiple times showing the dose information per group.

For more information about CTDIvol and DLP, refer to the Dosimetry section in the Technical Reference Manual for your system.

The Dose Report Screen Save can be networked, archived and filmed as part of the patient's record. The Dose Report can be set to be automatically auto filmed in the Film AutoFilm set up page on the View Edit screen.

\* Note that Enhanced Dose Reporting is not available on Single Processor systems due to hardware limitations. This feature is available for customer purchase by upgrading the CT operating console to Xstream GOC3 or GOC4 running Q4 2005 Common Software Platform.

**Enhanced Rotor Management:**

This feature optimizes rotor usage. The change in rotor management will help prolong X-ray tube life.

This feature is enabled automatically if the system recognizes that a GE supplied X-ray tube is installed. This feature is a purchased option for tubes not supplied by GE.

**Data Logging:**

ProDiags Enhancement to include tube usage and RTS. (See section 4 Data Logging Setup)


**Quality fixes included:**

The following quality problems are fixed in patch G.6:

- Reported DLP (dose length product) at the end of the exam is incorrect when returning to view edit to review dose information when AutomA is used.
- Failure of second scout in the next exam while recon from prior exam is still in progress.
- Fixes for service tools for investigating DAS issues.
- Fix for clean up of logs preventing them from growing too large.



### ***FMI Install Process Flow***

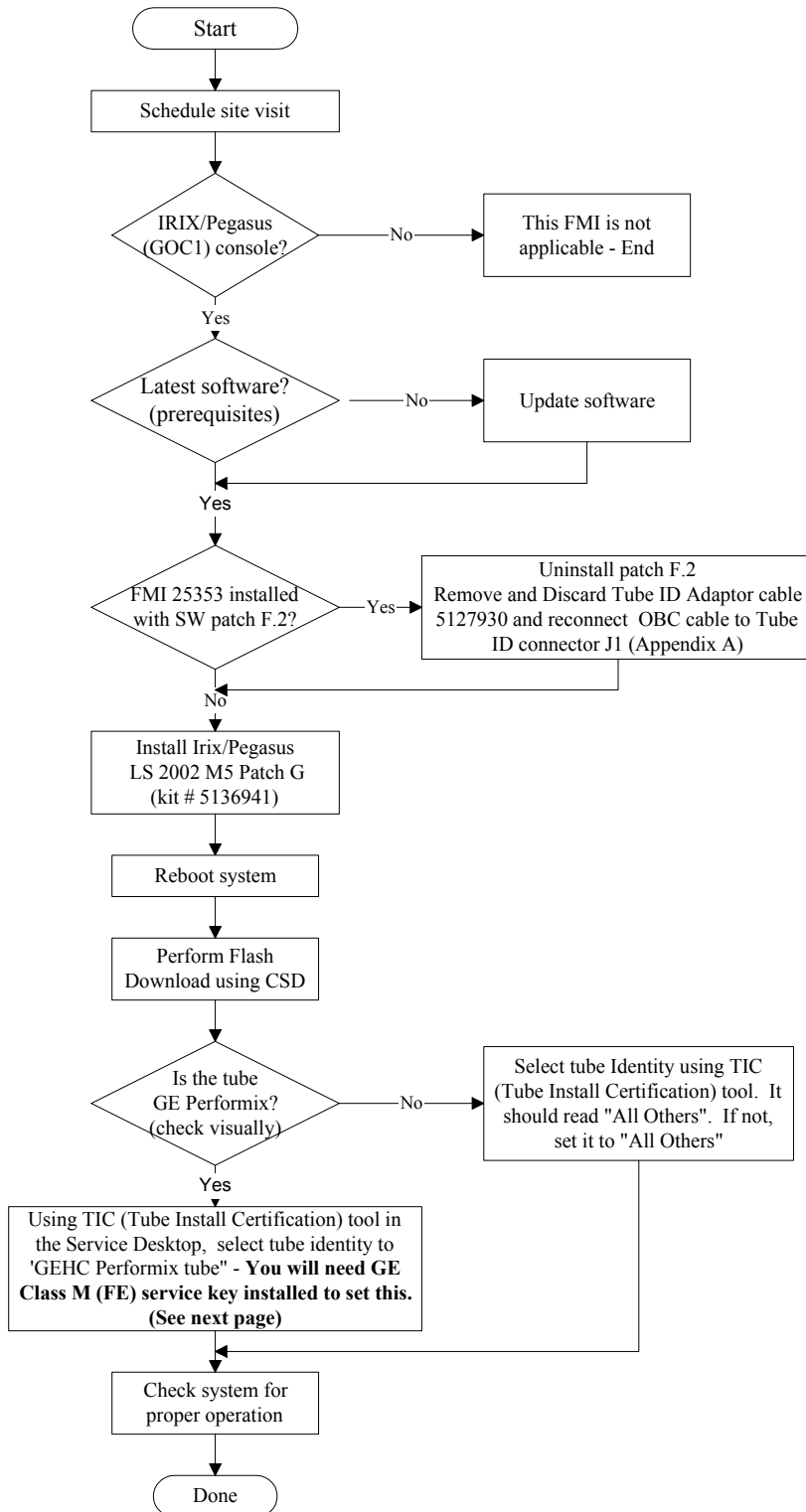



# Procedure

## 1 Preliminary Checks

1. Check to make sure that the system has the latest software. The latest software versions are listed in the Prerequisites section of this document. If the system does not have the latest software, please update the software before installing this FMI

### Notice:

- ***Install this FMI only in Irix/Pegasus (GOC1) systems.***
- ***Do not install this FMI to CT/PET systems***

2. Verify and record patches currently installed.

Open Unix Shell

Type **more /usr/g/patch/inst\_patches** <Enter>

 (Record patch info here.)

3. **If the system has been updated with software patch LS2002 M5 Patch F.2, Stop here. Execute instructions given in Appendix A to remove Patch F.2 before installing Patch G.6.**

## 2 Installation of Software Patch LS2002 M5 Patch G.6

Installation of Patch G.6: You need the CDROM LS2002 M5 Patch G.6

1. Shut down all CT Applications
2. Open a Unix Shell.
3. Type **su -** <Enter>
4. Type password **#bigguy** <Enter>
5. Insert LS2002 M5 Patch G.6 CD into the CT CD drive. ( CD # 5142334)
6. Type **mountCD /mnt** <Enter>
7. Type **cd /mnt** <Enter>
8. Type **./install\_patch -p ls2002\_m5\_patch\_G\_6** <Enter>

Patch will start loading, information will be displayed regarding the patch loading.  
This process takes ~ 2 minutes.  
Text will end with "POST-INSTALL"

9. Type **cd /** <Enter>
10. Type **umount /mnt** <Enter>
11. Type **eject 1 6** <Enter>

12. Remove the CD from the disk drive
13. Type **reboot** <Enter>  
System will go down and then come back to applications.
14. When back to Applications level, open a UNIX shell
15. As ctuser type **index\_img\_uid** <Enter> (This step is indexing the image table).
16. Confirm that the patch has been successfully installed by typing:

**showprods | grep ls2002\_M5\_patch\_G\_6** <Enter>

The output will be similar to:

*/ ls2002\_M5\_patch\_G.6 02/05/2005 Patch G.6 for ls2002 M5*

*/ ls2002\_M5\_patch\_G.6.img 02/05/2005 Patch files for ls2002 m5*

*/ ls2002\_M5\_patch\_G.6.img.ocroot 02/05/2005 Files that go onto ocroot.*

17. Close Unix Shell
18. Open Common Service Desk Top (CSD) and run Flash Download.  
Select Utilities  
Select Install Folder  
Select Flash Down Load Tool  
Select Update  
Dismiss when done
19. Select Pink Shutdown button  
Select Restart
20. When system returns to Applications, select tube identity following step 3.

#### PLEASE NOTE

- After each re-boot, 'un-recognized X-Ray Tube' message will be displayed as shown in the Appendix C since the tube identity has not been selected. Also displayed on CSD Homepage.
- If there is a need to uninstall software G.6, see Appendix B for instructions.


### 3 Selection of Tube Identity

Select Tube Identity ('GE Medical Systems Performix Tube' or 'All Others') using the Tube Install Certification Tool. The tool is in Common Service Desk Top under **Configuration**. Follow the process shown in the Process Flow on the following page.

**Verification of GE Medical System Performix tube: This can be verified by checking the tube insert model number (2120785 or 2120785-2), which can be found on tube rating plate.**

**Note:** Do not double click the Tube Certification Tool. This will cause more than one tool window to be displayed, one with message 'service key information not available'. If this happens, close window and try again. If this happens more than once, re-boot system and try again.

## Selection of Tube Identity

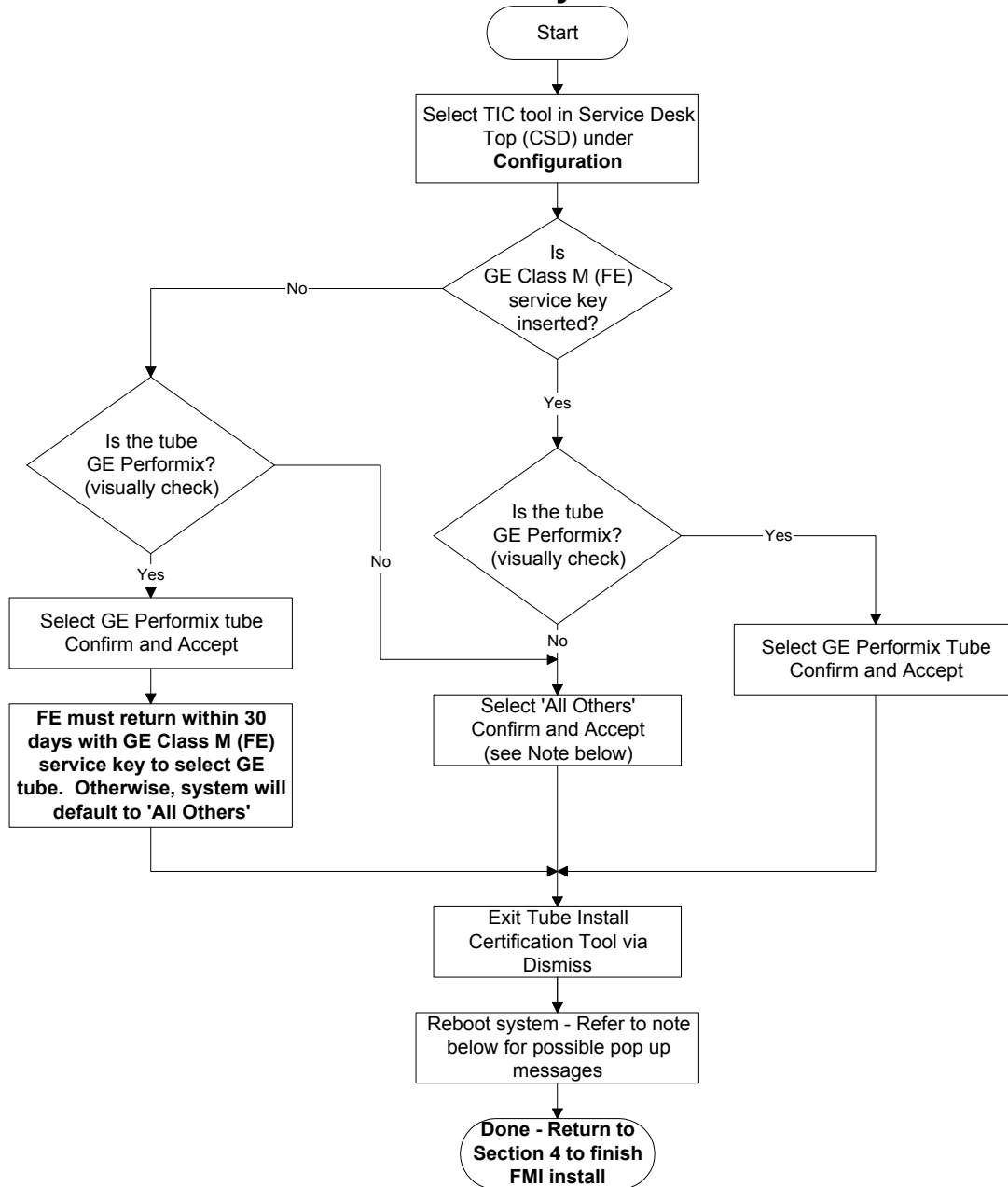

Note: Possible Pop Up messages at start up:

1. If the tube selection is GE tube without time limit, no pop up messages relating to tube identity on boot up.
2. If the tube selection is GE with time limit, a message with days remaining for GE FE to verify tube will be shown on boot up until GE FE removes the time limit with service key. If not verified by GE FE within the time limit, the system will default to 'All Others'.
3. If the tube is set to 'All Others', un-recognized x-ray tube installed message will always be shown at boot up.

## 4 Data Logging Setup:

Patch G.6 includes data logging enhancements that will push data to headquarters for analysis. These tasks are enabled for Octane Dual Processor systems if Class M software is installed and connectivity exists. These new tasks are PushTubeUsage and PushRTS.

This push of data is NOT enabled for Single Processor systems and should NOT be scheduled for Single Processors due to processor limitations and the affects on system performance.

Additionally, for both Dual and Single Processor systems, two existing tasks, (idm & idm\_send\_data) require schedule changes.

### 1. To check processor configuration:

1. Open a Unix Shell.
2. Type **check\_config** <Enter>

**Printout will look similar to the following:**

```
Product: LightSpeed16
sw_version: 405I.2_H4.0M5
machine_type: Octane
host_id: 56037637
tube_type: 15
gantry_type: GANTRY_COMPACT
No. of processors: 2 (This line will indicate processor configuration, Single =1 or Dual = 2).
Processor speed: 400
memory: 1536
```

Step-by-Step instructions for ProDiags adjustments are provided here, for a screen-shot of ProDiags screen, see page11.

### 2. ProDiags Adjustments (Single and Dual Processor Systems)

1. Open a Unix Shell.
2. Type **prodiags** < Enter >
3. Select **Schedule** Tab
4. Go to Task **idm**, change Hour to **4 PM** and change Minute to **35**.
5. Go to Task **idm\_send\_data**, change Hour to **4 PM** and change Minute to **40**.
6. Click **Apply** at the bottom of the screen.
7. View Tasks **PushTubeUsage** and **PushRTS** at bottom of schedule these tasks should be listed and scheduled for Dual Processor systems and should NOT be listed for Single Processor Systems. Please see Note below concerning changing schedule for these tasks.

**NOTE:** The default schedules for these tasks provided in this release are optimized for most sites are highly recommended. Careful planning is required before making any schedule changes as the amount of data pushed can vary based on the frequency. Changing the schedule may result in undesirable changes in system performance and/or inadequate data collection.

8. Close Proactive Diagnostics window with pull-down in upper left hand corner.

If Software is reloaded on the system at any time, tasks idm and idm\_send\_data should be re-adjusted to the times set above.

### ProDiags Screen for Octane Dual Processor Systems:

**Proactive Diagnostics**

Schedule Task Time Slot Log & Result

**Schedule Management**

Schedule Management provides the following functionalities

1. Add a task to the schedule
2. Delete a task from the schedule
3. View the schedule

| Task          | Iterations | Day       | Hour   | Minute | Type       | Time Slot |
|---------------|------------|-----------|--------|--------|------------|-----------|
| PassChg       | 1          | Daily     | 05 AM  | 30     | Background |           |
| pd_HouseK...  | 1          | Daily     | 04 AM  | 30     | Background |           |
| pd_ASC_No...  | 1          | Daily     | Hourly | 28     | Background |           |
| healthpg      | 1          | Daily     | 04 AM  | 45     | Background |           |
| IIP_HouseK... | 1          | Daily     | 06 AM  | 12     | Background |           |
| IIP_HouseK... | 1          | Daily     | 12 PM  | 12     | Background |           |
| IIP_HouseK... | 1          | Daily     | 06 PM  | 12     | Background |           |
| IIP_HouseK... | 1          | Daily     | 12 AM  | 12     | Background |           |
| idm_send...   | 1          | Daily     | 07 AM  | 00     | Background |           |
| idm           | 1          | Daily     | Hourly | 00     | Background |           |
| ctrdi         | 1          | Daily     | 12 PM  | 59     | Background |           |
| PushTubeU...  | 1          | Wednesday | 10 AM  | 16     | Background |           |
| PushRTS       | 1          | Daily     | 10 AM  | 07     | Background |           |

New Delete Reload Apply Help

Status: Reading data... Done

**NOTE: If sw is reloaded in future, the schedule for these tasks needs to be reset!**

Set Time for **idm\_send\_data** task to: **16:40 or 4:40PM**

Set Time for **idm** task to : **16:35 or 4:35 PM**

**NOTE: New Tasks added to push data to backoffice for Dual Processor Systems.**

## 5 Enhanced Dose Reporting

No action is required on Patch installation, the following material is provided for informational use only.

Systems with a **dual processor** configuration will have an Enhanced Dose Reporting feature. The Dose reporting tool will create a screen save image for the dose information saved in series 999. This image can be filmed, networked and archived.

For systems with a **single processor** configuration, the screen save image is not available.

## 6 System Checkout

Perform the following scans (Abdominal or Head)

50 images at 0.5 sec rotation at 120KV, 40mA

The system should scan and image without any errors

## 7 Finishing Up

1. Make sure that the 'Tips and Workaround' document, 5115418-100 supplied with this FMI kit is left with the customer. Also, leave the Customer Letter, 5142333-147 with the customer.
2. Close FMI per regional procedure.
3. Remove this FMI document from the customer site.

## Appendix A

### Removal Of Patch F.2

#### Patch F.2 Uninstallation:


Note:

If you uninstall patch F.2, you should also remove cable adaptor, 5127930 as shown following the uninstall instructions

1. Shutdown all CT applications
2. Insert the LS2002 M5 Patch G.6 into the CD drive
3. Open a unix shell
4. Type **su** - <Enter>
5. Type password **#bigguy**
6. Type **mountCD /mnt** <Enter>
7. Type **cd /mnt** <Enter>
8. Type **./uninstall\_patch** <Enter>  
This will display a menu with all the patches installed and will ask you to enter the number for the patch that you want to uninstall.
9. To uninstall patch F.2, **enter the number specified for F\_2** <Enter>
10. Type **cd /** <Enter>
11. Type **umount /mnt** <Enter>
12. Type **eject 1 6** <Enter>
13. Remove CD out of the CD drive
14. Shutdown system and power off
15. Remove cable adaptor, 5127930 using the procedure outlined in the following pages
16. After cable is removed, bring system back up to applications
17. Once the system has rebooted run Flash Download (see step 15 of the install Instructions)
18. Select Pink Shutdown button  
Select Restart

## Appendix A – Contd.

### Removal of Cable Adaptor, 5127930

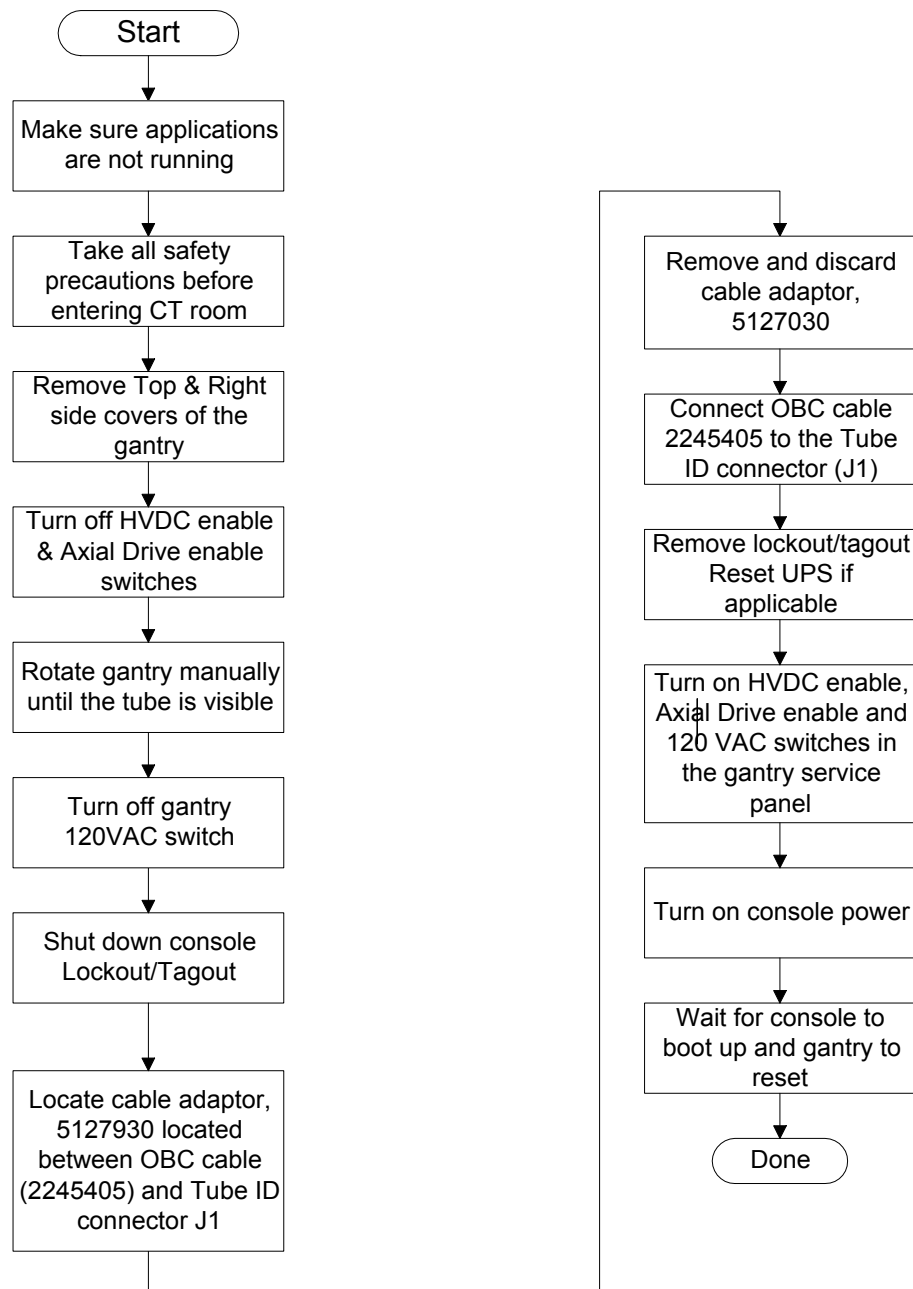



## Appendix A – Contd.


### Removal of Cable Adaptor, 5127930

1. Make sure that applications are not running
2. Take all safety precautions before entering the scan room
3. Lower table to home (lowest) position
4. Remove gantry top and side covers

**WARNING**

FAILURE TO DISABLE THE AXIAL DRIVE MAY RESULT IN SERIOUS INJURY DUE TO UNEXPECTED GANTRY ROTATION.

5. Turn off high voltage DC enable switch (HVDC) and the Axial Drive enable (AX DR) switch to off position in the service panel inside the gantry

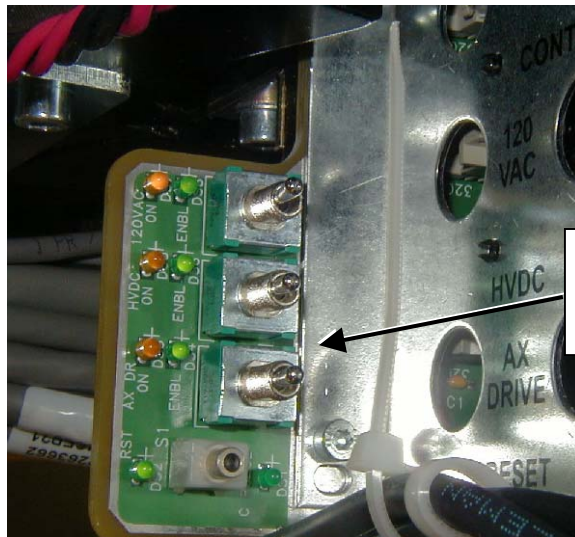


HVDC enable,  
Axial Drive enable and  
120VAC switches

6. Rotate the gantry with care until the tube is visible.
7. Now, turn off the 120VAC switch in the service panel
8. Shut off power to the console and follow lockout/tag out procedure to make sure that power is off while making the changes to the gantry
9. Locate cable adaptor, 5127930 which is connected to OBC cable, 22445405 at one end and to the Tube ID connector J1 at the other end. It is behind the tube cooling tank power cable.
10. Remove the cable adaptor, 5127930 and discard.
11. Connect cable 2245405 to the tube at connector J1.
12. For easy access to tube connector J1, if you have to disconnect Fan & Pump power cable, **REMEMBER TO RECONNECT AT COMPLETION.**

## Appendix A – Contd.

### Removal of Cable Adaptor, 5127930

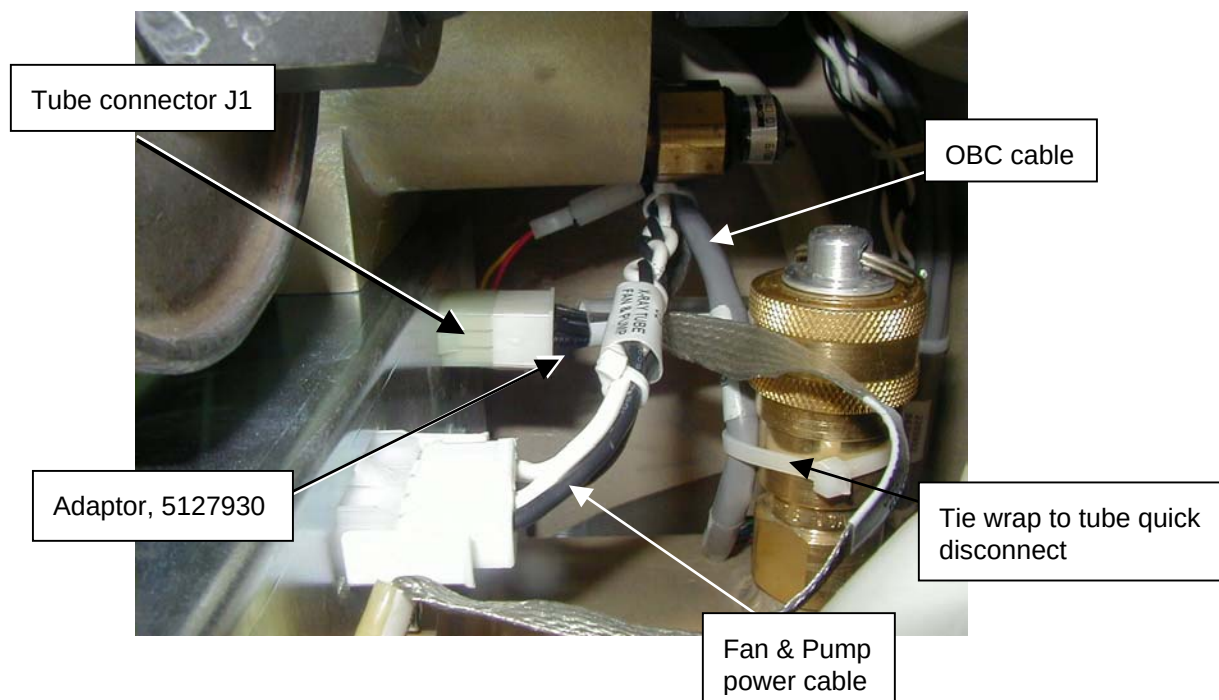

IMPORTANT: Visually check the gantry to make sure that:

- ***The power cable to the tube cooling tank is connected***
- ***The tube ground wire is in place***

13. Remove lockout/tag out so that power can be applied to the system. If the site uses UPS, reset UPS.
14. Turn on all three switches (Axial Drive, HVDC and 120VAC) in the gantry
15. Power up console
16. Wait for the console to boot up and the gantry to reset
17. Return to section 2, Installation of Software Patch LS2002 M5 Patch G.6

## Appendix B


### Patch G.6 Uninstallation:

Note:

If you uninstall patch G.6, you need to restore system to the original system configuration. Patches, which were installed before G.6 should be restored.


1. Shutdown all CT applications
2. Insert the LS2002 M5 Patch G.6 CD into the CD drive
3. Open a unix shell
4. Type **su** - <Enter>
5. Type password **#bigguy**
6. Type **mountCD /mnt** <Enter>
7. Type **cd /mnt** <Enter>
8. Type **./uninstall\_patch** <Enter>  
This will display a menu with all the patches installed and will ask you to enter the number for the patch that you want to uninstall.
9. To uninstall patch G.6, **enter the number specified for G\_6** <Enter>
10. Type **cd /** <Enter>
11. Type **umount /mnt** <Enter>
12. Type **eject 1 6** <Enter>
13. Remove CD out of the CD drive
14. Type **reboot** <Enter>
15. Once the system has rebooted run Flash Download (see step 15 of the install instructions)
16. Select Pink Shutdown button
  - i. Select Restart

## Appendix C - Screen Displays

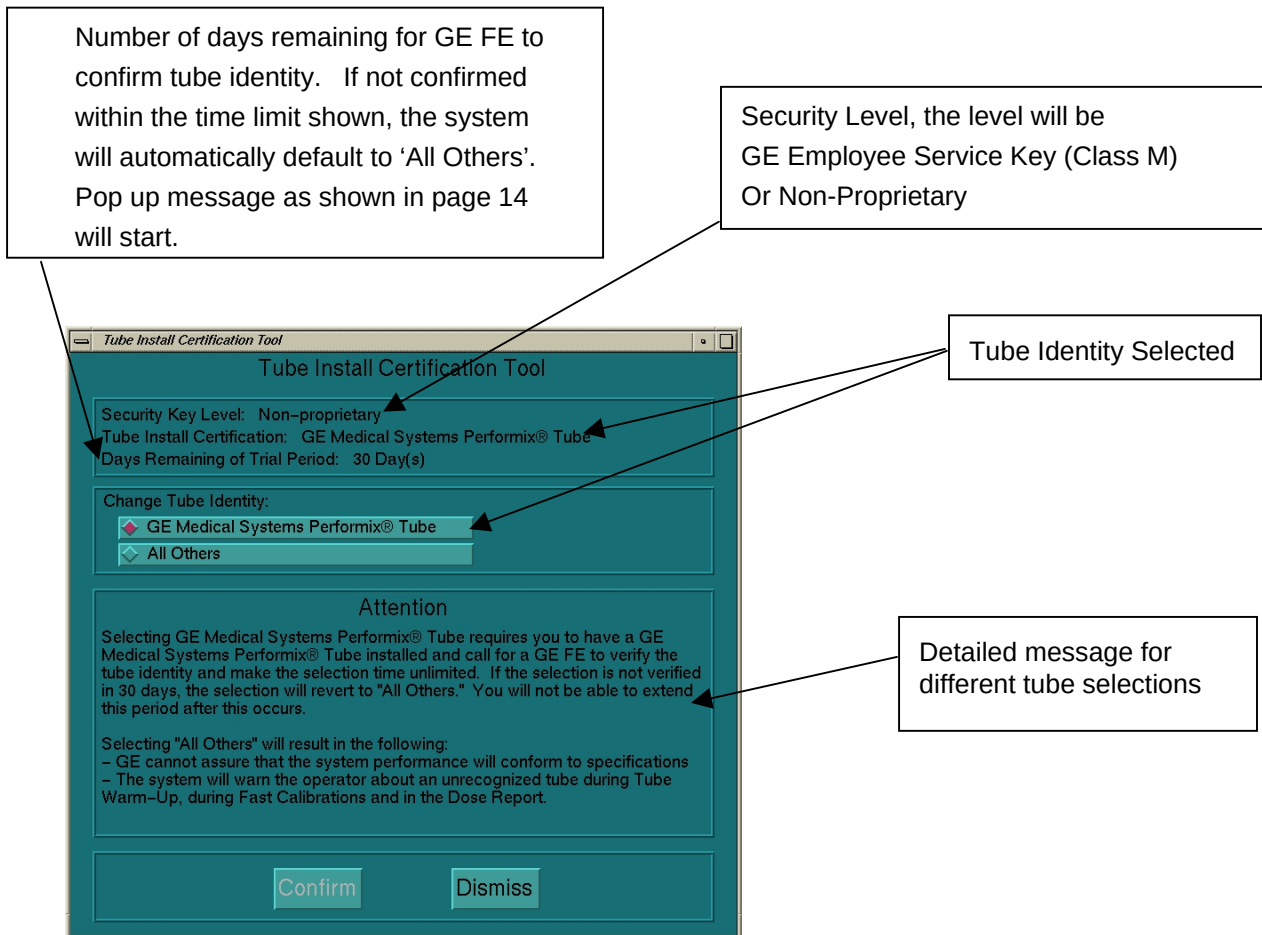


Example of Tube Install Certification Tool

## Appendix C - Screen Displays (Contd.)

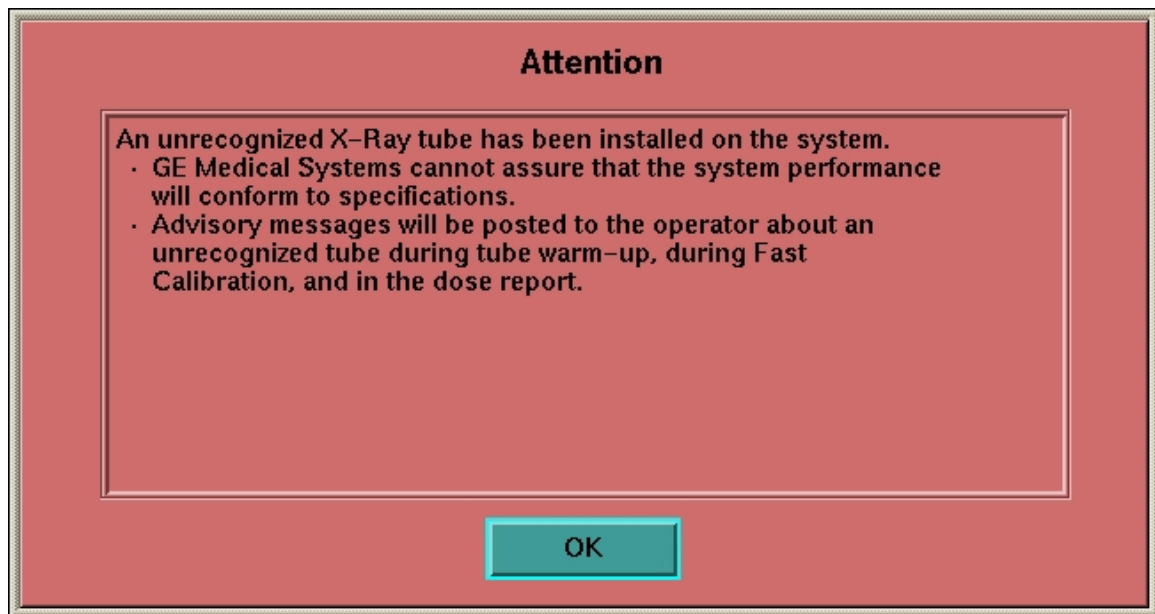

**Message at each boot up when unrecognized X-Ray Tube is installed in System**

## Appendix C - Screen Displays (Contd.)

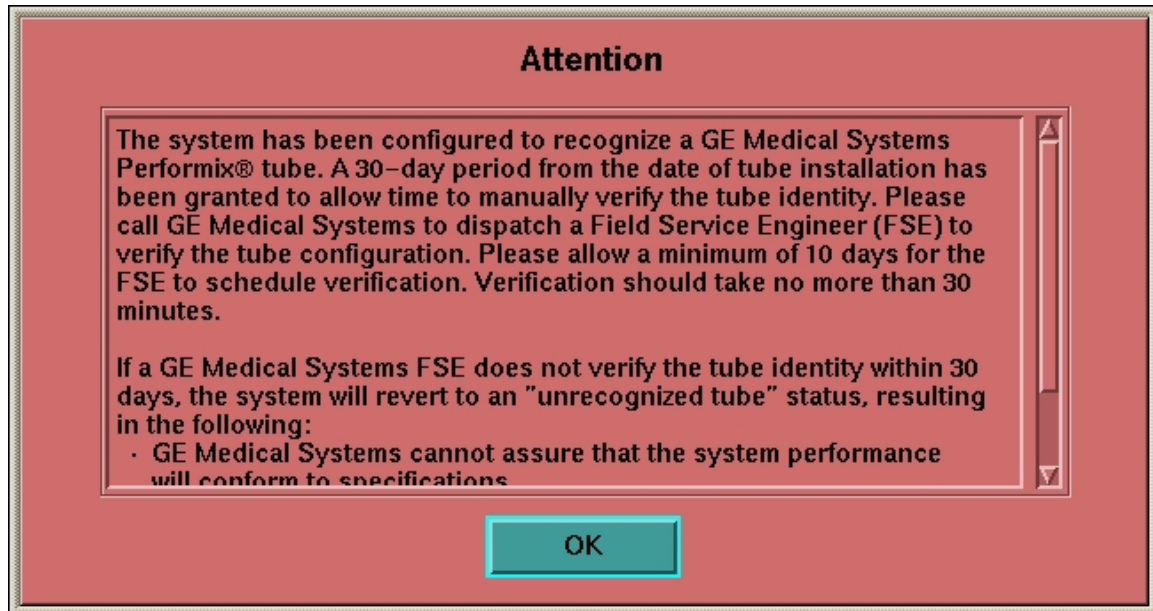

**Message at each boot up when GE Tube is selected but not verified by  
GE Field Engineer**

## Appendix C - Screen Displays (Contd.)

|                         |         |                    |                  |                 |               |
|-------------------------|---------|--------------------|------------------|-----------------|---------------|
| Patient Name: Sam Ample |         |                    | Exam no: 507     |                 |               |
| Accession Number:       |         |                    | 21 Jun 2005      |                 |               |
| Patient ID: Sample      |         |                    | LightSpeed Plus  |                 |               |
| Exam Description:       |         |                    |                  |                 |               |
| Dose Report             |         |                    |                  |                 |               |
| Series                  | Type    | Scan Range<br>(mm) | CTDIvol<br>(mGy) | DLP<br>(mGy-cm) | Phantom<br>cm |
| 1                       | Scout   | -                  | -                | -               | -             |
| 2                       | Helical | S30.000-S7.500     | 29.75            | 106.93          | Body 32       |
| 2                       | Helical | I72.500-I95.000    | 29.75            | 106.93          | Body 32       |
| Total Exam DLP:         |         |                    |                  | 213.86          |               |
| 1/1                     |         |                    |                  |                 |               |
| W:1 L:-1025             |         |                    |                  |                 |               |

Dose Report Screen when GE Medical Systems  
X-Ray Tube is installed in System

## Appendix C - Screen Displays (Contd.)

|                                                                                                                                                                                                                                                                                                                                |         |                    |                  |                 |               |
|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------|--------------------|------------------|-----------------|---------------|
| Patient Name: Test Name                                                                                                                                                                                                                                                                                                        |         |                    | Exam no: 504     |                 |               |
| Accession Number:                                                                                                                                                                                                                                                                                                              |         |                    | 21 Jun 2005      |                 |               |
| Patient ID: sfsf                                                                                                                                                                                                                                                                                                               |         |                    | LightSpeed Plus  |                 |               |
| Exam Description:                                                                                                                                                                                                                                                                                                              |         |                    |                  |                 |               |
| Dose Report                                                                                                                                                                                                                                                                                                                    |         |                    |                  |                 |               |
| Series                                                                                                                                                                                                                                                                                                                         | Type    | Scan Range<br>(mm) | CTDIvol<br>(mGy) | DLP<br>(mGy-cm) | Phantom<br>cm |
| 1                                                                                                                                                                                                                                                                                                                              | Scout   | -                  | -                | -               | -             |
| 2                                                                                                                                                                                                                                                                                                                              | Helical | S30.000-S7.500     | 29.75            | 106.93          | Body 32       |
| Total Exam DLP:                                                                                                                                                                                                                                                                                                                |         |                    |                  | 106.93          |               |
| Unrecognized tube in use - Dose may vary. The reported dose information is calculated based on empirical observations of systems with GE Medical Systems Performix® tubes. GE cannot assure the accuracy of reported dose information for any configurations that include tubes other than GE Medical Systems Performix tubes. |         |                    |                  |                 |               |
| 1/1                                                                                                                                                                                                                                                                                                                            |         |                    |                  |                 |               |
| W:1 L:-1025                                                                                                                                                                                                                                                                                                                    |         |                    |                  |                 |               |

Dose Report Screen when un-recognized  
X-Ray Tube is installed in System