

**ADMINISTRATIVE INFORMATION****FMI No: 25365****FMI Title:** 06MW03.5 SW for CT- LS 1X/2X/3X/4X/5X & HS QX/i with GRE Consoles

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS – AM:	Issue Date:
☐ Mandatory	GEMS – E:	
☐ Mandatory on Request	GEMS – A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS – Americas	GEMS - Europe	GEMS – Asia
GE Medical Systems On-Line Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:

Note: Critical tube usage information is converted by 06MW03.5 LFC procedure. **It is important to execute the LFC instruction set as provided to ensure all warranty information is properly converted and GE Healthcare's financial obligations are maintained.**

Perform software Load using the procedure outlined in service CD-ROM 5167277-200.

- Open CD-ROM 5167277-200, *FX/Xtream Workflow w/Q1 '06 SW - Upgrade Procedures*
- Click on continue
- Click on the link to LFC Procedure – Version 06MW03.5 on the right pane
- Follow the instructions for software installation

Unit Tracked:**FMI Instructions Part No.: 5161945-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas

GEMS – Europe

G.E. Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe
FMI Admin. 322
BP34
F78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25365

FMI Title: 06MW03.5 SW Update for CT- LS 1X/2X/3X/4X/5X & HS
QX/i with GRE Consoles

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments

Effectivity

The FMI effectivity includes CT GRE/Xtream (GOC3 and GOC4) console systems for:

- LightSpeed 1X
- LightSpeed 2X
- LightSpeed 3X
- LightSpeed 4X
- LightSpeed 5.X
- HiSpeed QX/i

Model No.

- 2377708
- 2377708-4
- 2377708-7
- 2377708-10
- 2377708-13
- 2377708-20
- 2377708-21
- 2377708-22
- 2377708-30
- 2377708-31
- 2377708-32
- 2377708-40
- 2377708-41
- 2377708-42
- 2377708-43
- 2377708-50
- 2377708-51
- 2377708-52
- 2377708-53
- 5111750

Purpose

The purpose of this FMI is to upgrade system with the new 06MW03.5, software on LightSpeed 1X/2X/3X/4X/5X and HiSpeed QX/i systems with GRE (GOC3 and GOC4) consoles.

The software fixes several key quality and safety issues. Also introduces several new features.

Prerequisites

Software:

Please make sure that the software revisions are latest for your system (as shown below) before loading the new software.

SYSTEM	SOFTWARE VERSION
LightSpeed 1X	04MW13.7; 04M44.X; 05MW14.5
LightSpeed 2X	04MW13.7; 04M44.X; 05MW14.5
LightSpeed 3X	308L.5_H3.3M4; 308L.8_H3.3M4; 04MW13.7; 04MW44.X; 05MW14.5
LightSpeed 4X	405L.5_H4.2M4; 405L.8_H4.2M4; 04MW13.7; 04MW44.X; 05MW14.5
LightSpeed 5X (HP100, HP80 and RT)	04MW13.7; 04MW44.X; 05MW14.5
LightSpeed 5X (HP60)	04MW23.1; 04MW44.X; 05MW14.5
HiSpeed QXi	103L.4_HSQXIM4; 04MW13.6; 04MW44.X; 05MW14.5

Note:

Software versions 405L.8 and 308L.8 were deployed by the last globally distributed FMI's 25338 and 25337 respectively.

Furnished Materials

Below are the items included in kit# 5137245.

Item #	Qty	Description	Part #
1	1	FMI Document	5161945-100
2	1	Software Collector with rating plate	5167609
3	1	HP60 AX Motor HV Sticker	5146059
4	2	Terminal 20 AMPS 600 Volts	2275176
5	1	<i>FX/Xtream Workflow w/Q1 '06 SW - Upgrade Procedures</i> CDROM (includes 06MW03.5 LFC)	5167277-200
6	2	DVD RAM Media, 9.4 GB (Blank Disk)	2337414
7	1	Advanced Applications on OC OP Info	2404359-247
8	1	InSite™ Update (Software Download) Marketing Kit	5164884
9	1	Tube Usage Converter Tool for SW Rls 06MW03.5	5172979
10	1	New Features Guide for 04MW13.7 and 04MW44.9 Software. Tips and Workaround	5166132-100
11	1	CT Light Speed series CD	2369740-247
12	1	Customer Letter	5142333-147
13	1	Common Platform Q1 2006 Release Applications Tips and Precautions - English	5167065-100

Tools Required

LightSpeed Class C Software (tab 182 distribution)

Part # 5165729

LightSpeed Class M Software (tab 182 distribution)

Part # 5154540

Standard Service Tool Kit

Customer purchased options

Personnel Requirements

Travel: 1 Hour

Labor: 6 Hours

Field Engineer: 1

Product Safety Fixes

- Scout Image Fails to recon

This was caused by characters added at the end of the Exam Description when the shift or control key was pressed along with the enter key.

- Additional slices added when scans are prescribed with gantry tilt

Improvements have been made in Graphic Rx and ScanRx

- Show Localizer lines are in wrong location if user skips series quickly

- Smart Step Images out of sequence

- Auto mA changes automatically or reduces to 10 mA

This was due to how Noise index was updated when the user changed slice thickness. Prior to this patch, the Noise index adjustments were made from the Reference Noise Index. The change is to adjust Noise Index on Slice Thickness change with the Rx'd Noise Index. This will assure that when the user changes slice thickness and then returns to the original slice thickness the Noise Index will default to the original value.


- System recommends optimize to 35 mA
- Show Localizer lines are in wrong location if user skips series quickly
- DLP is higher than expected caused by selecting Repeat Series Quickly.

The dose reported was incorrect.

Improvements

- Apostrophe in Hospital Name causes Tube Warm Up to fail
- SmartPrep Diagnostic delay not as expected.

The diagnostic delay is adjusted at Confirm to account for table movement required between the Monitor and Scan phase locations. If the Diagnostic Delay prescribed is less than what the system can do, the expected delay will be posted at Confirm.

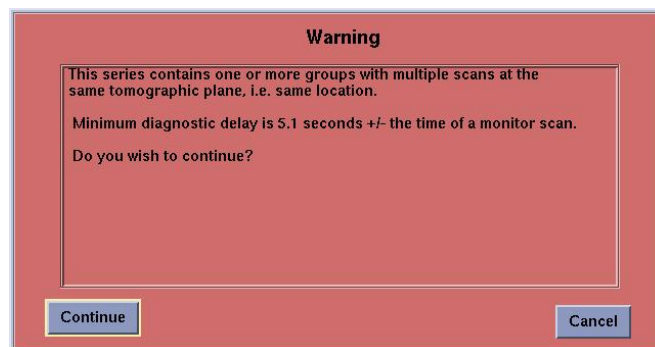

- SmartPrep fails to start
- Improvements have been made in communication between the ScanRx and Display Processes to assure reliable use of SmartPrep
- Scan aborts due to Axial not ready

- Start Locations have odd numbers.
Adjustments to Start Locations from Graphic Rx will be rounded off in .25mm increments. This will give more meaningful scan locations for the user.
- Display hung with 10% loading
- System Hung with Can not dismiss Attention Window
- Images can't be retrieved from PACS
- System shuts down with a ScanRx crash
- Dcserver crashes with Storage Commitment requests
- Auto Voice fails to play if the table is moved during a long prep delay
- Images are filmed cut off
This was caused when filming images with GSE from a non-square viewport.
- Update Apache HTTP server to be more secure
- Clean up Pro Diags so Reli Monitor does not fail
- Eliminates scan aborts caused by Axial Motor control logic
- Eliminates scan aborts caused by Rotor not at speed on HP60 systems.

AVA (Advanced Vessel Analysis) Measurements Wrong after Edit Contour

AVA measurements wrong after edit contour – The Dmin and Dmax values were recalculated incorrectly after editing the contour. The values are now calculated correctly.

When using AVA, the user should be aware of the Measurement accuracy of the tool. Please note these new accuracy guidelines when performing measurements after using the manual edit contour tool.

Measurement Accuracy

CAUTION: Accuracy of cross-section distance measurements should be +/- 10% for vessels with a diameter between 10 and 20 voxels. It should be +/- 5% for vessels with a diameter greater than 20 voxels.

CAUTION: Accuracy of cross-section area measurements should be +/- 40% for vessels with a diameter between 10 and 20 voxels. It should be +/- 20% for vessels with a diameter greater than 20 voxels.

CAUTION: Accuracy of volume measurements should be +/- 40% for vessels with a diameter and a length between 10 and 20 voxels. It should be +/- 20% for vessels with a diameter and a length greater than 20 voxels.

CAUTION: Accuracy should be better than +/- 10 degrees for angle measurements.

New Features (Software 06MW03.5):

SmarTube

There are further enhancements to the SmarTube feature in this release. Additional software has been added to make a selection for tube identification based on tube type,

1. GE-Supplied x-ray tube (Performix™) and
2. All other tube types.

If a Non-GE-Supplied tube is found, a warning dialog box will be posted at system start-up.

Advisory messages will be posted on systems with Non-GE-Supplied tubes during Fast Cal and on the Dose Report screen save image indicating that information may not be totally accurate because the full characteristics of the tube may not match those of a GE-Supplied x-ray tube. If the installed x-ray tube is a GE-Supplied tube but the field service engineer did not have all the appropriate information to complete verification, a message will be displayed for 30 days or until the field service engineer is able to complete tube installation verification within the 30-day window.

If the verification is not completed in the 30 day window, software will revert to unrecognized tube state and messages for non-validated x-ray tube will be seen during FastCal and on the Dose Report screen save image.

Dose Report Screen Save:

Dose Report Screen Save feature has been implemented in this release. The Dose Report is a screen save image that is automatically saved to Series 999 at normal End Exam. The Dose Report is a copy of the information displayed in the Dose Information Area on the View Edit screen. The report includes the CTDIvol, Dose Length Product (DLP) and total DLP for the exam. Scouts will be listed, however, there is no recognized method for measuring dose for scouts so the information for these series will always be blank. If there is more than one group in a series, the series will be listed multiple times showing the dose information per group.

For more information about CTDIvol and DLP, refer to the Dosimetry section in the Technical Reference Manual for your system.

The Dose Report Screen Save can be networked, archived and filmed as part of the patient's record. The Dose Report can be set to be automatically auto filmed in the Film AutoFilm set up page on the View Edit screen.

Enhanced Rotor Management:

This feature optimizes rotor usage. The change in rotor management will help to prolong the x-ray tube life.

This feature is enabled automatically if the system recognized that a GE supplied x-ray tube is installed. This feature is a purchased option for tubes not supplied by GE.

Data Logging:

Enhancement to ProDiags is done to include tube usage and RTS.

InSite™ Update

The InSite™ Update (Software Download/Update Management) feature allows the customer to get software updates and small features via a broadband automatic update process. The system must be InSite™ broadband connected to maximize the full benefit of the automatic InSite™ Update feature. This tool is accessible from the Common Service Desktop (CSD).


For a detailed description and operation of this feature, please see the 'InSite™ Upgrade' section of the document, 5169826-200, Training Material, New Features Deployed by SW Release 06MW03.5 for VCT M4 and PF1.5 M3.

Product Network Filter (PNF)

Networked medical imaging systems such as CT can be vulnerable to security risks. This software version includes Product Network Filters (PNF), which addresses the concern of network security by providing a host-based firewall with the capability to filter and block unwanted network communication. The default configuration will be adequate for most all systems. For details regarding PNF and configuration information, please refer to the 'Product Network Filter' section of the document, 5169826-200, Training Material, New Features Deployed by SW Release 06MW03.5.

Refer to 5166132-100 for new features from the 04MW13.7, 04MW44.9 and 05MW14.5 software releases.

Process Flow

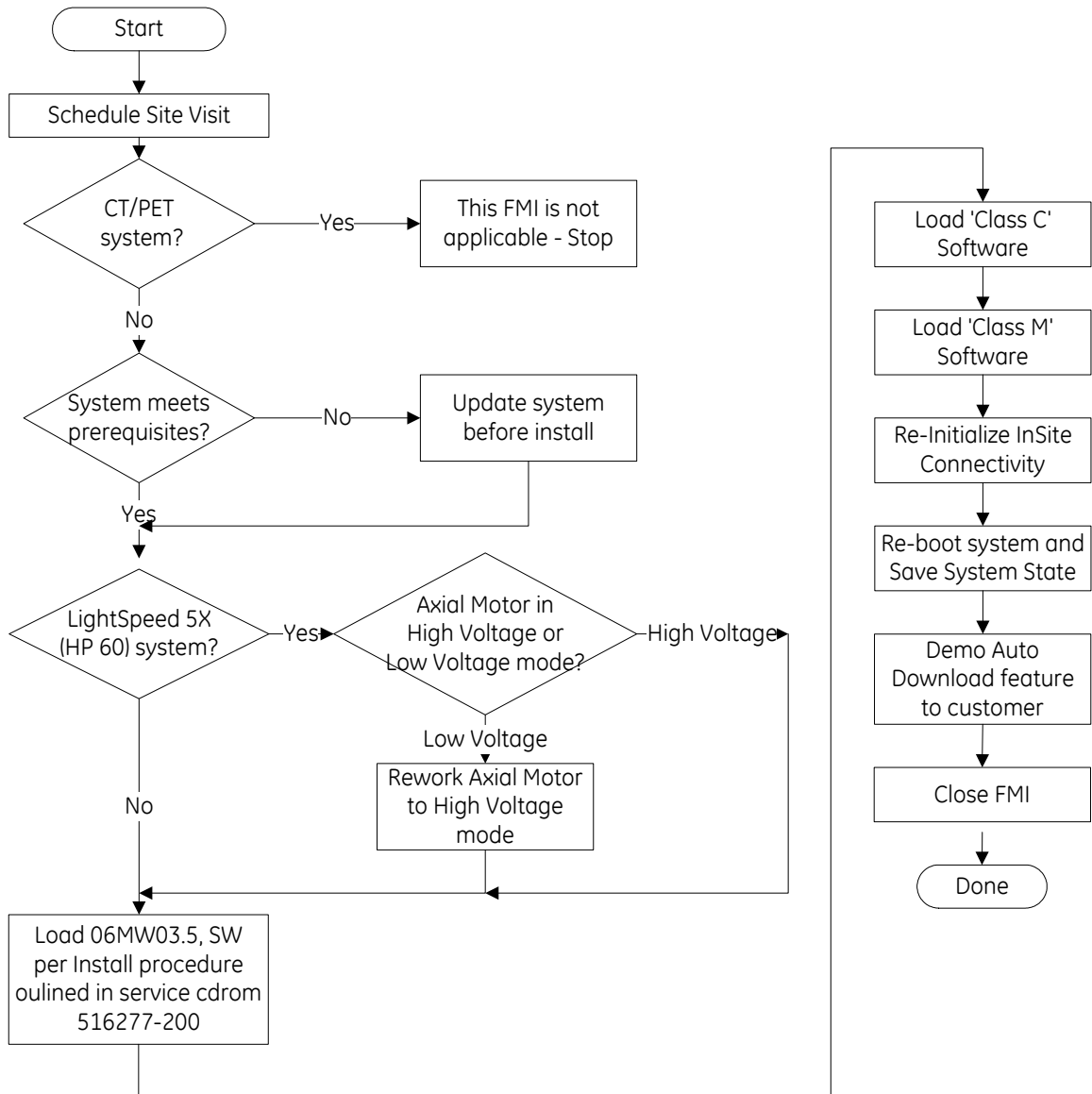

Procedure

1 Preliminary Checks

- 1 If this is CT/PET system, this FMI is not applicable. STOP.
- 2 Check the system's software version before executing this FMI.
Refer to the Prerequisites section at the beginning of this FMI instruction (page 2).
STOP if software prerequisites are not met. Load the proper software release before loading this FMI.
- 3 Open the Common Service Desktop Home tab, and record the Application software version listed under "Current System Status" here:

This information will be used later, in the Tube Usage Conversion Tool section.

Tube usage data dictates GEMS' warranty responsibilities to the customer.

2 Installed Options and InSite™ Connectivity Information

Before the upgrade, record the following information:

Installed options:

- Open Unix shell
- Type `swokinstall -p` <Enter>

From the swokinstall -p output, record system options installed in the system.

InSite™ Connectivity Information:

- Open Unix shell
- Type `cd /usr/g/insite` <Enter>
- Type `more .insiteINFO` <Enter>

Record connectivity configuration parameters.

3 Axial Motor Control Rework (ONLY FOR HP60 SYSTEMS)




WARNING

POTENTIAL FOR SHOCK.
VOLTAGE MAY BE PRESENT.
USE PROPER LOCKOUT/TAGOUT PROCEDURES WHEN WORKING IN THE
GANTRY.

- 3.1 If your system is an HP60 system, Axial Motor Control work **REWORK** may be needed. This modification is needed to fix error “260107214, gantry failed to stop” in certain conditions. To check if your system is an HP60 system:

Remove Gantry right side cover

Check the label on the support leg (see figure below)

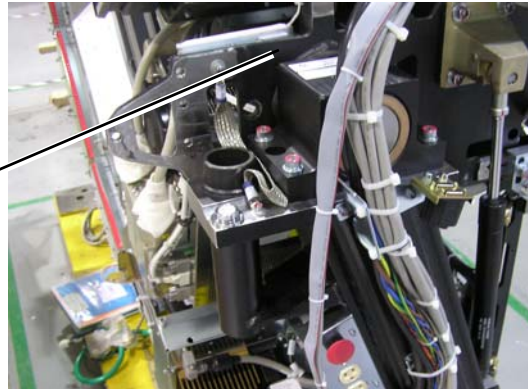

If your system is not HP60, skip this step and go to Section 4 ‘Software Load Procedure’.

If your system is HP60, check whether the rework has been done. To check if the HP system has the rework done,

Remove Gantry right Gantry side cover

Turn off the gantry service switches (Axial Drive and 120VAC)

Remove the support bar of the front cover

Check the side mesh cover of the Axial Motor. If there is a label saying ‘AMC H’, the rework has already been done. **Go to Section 4 ‘Software Load Procedure’.**

3.2 Modification of Wire Connections

Perform Axial Motor Control rework as outlined in the service CD-ROM 5167277-200:

- Open CD-ROM 5167277-200, Xstream FX 1Q '06 Upgrade CD-ROM
- Click on Pre-Install.
- Click on Axial Motor Control Rework.
- Follow the instructions for rework.

4 Software Load Procedure

NOTE:

Critical tube usage information is converted by 06MW03.5 LFC procedure. It is important to execute the LFC instruction set as provided to ensure all warranty information is properly converted and GEHC's financial obligations are maintained.

Perform software Load using the procedure outlined in service CD-ROM 5167277-200:

- Open CD-ROM 5167277-200, *FX/Xtream Workflow w/Q1 '06 SW - Upgrade Procedures*
- Click on Continue.
- Click on the link for LFC Procedure - Version 06MW03.5 on the right pane.
- Follow the instructions for software installation.

5 Load Class C and Class M Software

To ensure proper load of the InSite™ software, please follow the following sequence:

- Load C Software
- Load M software
- Re-initialize InSite™ Connectivity

If the Class C and Class M Software install log exhibits errors, usually a repeat attempt at loading the software succeeds.

5.1 Load Class C Software

Load Class C software as follows:

1. Open a shell.
2. Enter root as superuser:
 - a. At the ctuser prompt, type: `su -` <Enter>
 - b. At the Password prompt, type: `#bigguy` <Enter>
3. Insert the Class C software CD into HP computer DVD Drive. Wait for the disk read to STOP.
4. Install the Class C software on the system. At the root prompt, type: `install.classC` <Enter>
5. Close window after install.
6. Verify Class C is properly loaded on the system.
 - a. Open a Unix Shell
 - b. At the root prompt, type: `cd /var/adm` <Enter>
 - c. At the root prompt, type: `ls` <Enter>
 - d. Find the `install.classC.log` file. If file found, this means the Class C software load completed.
7. Remove the Class C software CD from HP computer DVD Drive.

5.2 Load Class M Software

If your system does not have InSite™ Entitlement, go to section 7 (on page 17).
Load provided Class M software as follows:

1. Open a shell.
2. Enter root as superuser:
 - a. At the ctuser prompt, type: `su -` <Enter>
 - b. At the Password prompt, type: `#bigguy` <Enter>
3. Insert the Class M software CD into the HP computer DVD drive. Wait for the disk read to STOP.
4. Install the Class M software on the system. At the root prompt, type: `install.classM` <Enter>
5. Close window after install
6. Verify Class M is properly loaded on the system.
 - a. Open a Unix Shell
 - b. At the root prompt, type: `cd /var/adm` <Enter>
 - c. At the root prompt, type: `ls` <Enter>
 - d. Find the `install.classM.log` file. If file found, this means the Class M software load completed.
7. Remove the Class M software CD from HP computer DVD Drive.

6 Re-Initializing InSite™ Connectivity

6.1 Restore System State To Retrieve Class M Configuration Parameters

Restore 'Reconfig' part of the Restore System State to recover all previous InSite™ and Class M configuration parameters. Be sure to use the new State DVD with converted tube usage files.

6.2 Enable InSite™ Interactive Platform Configuration Utility

Enable the InSite™ Interactive Platform Configuration utility as follows:

1. Open a shell.
2. Enter root as superuser:
 - a. At the ctuser prompt, type: `su -` <Enter>
 - b. At the Password prompt, type: `#bigguy` <Enter>
3. Launch the IIP Admin Configuration utility. At the root prompt, type: `iipadmin config` <Enter>

The system displays the InSite™ Interactive Platform Configuration utility with the Proprietary Agreement Acceptance Screen.

4. Click on ACCEPT three times to display the Information Screen.
5. Click on the **ProDiags** Tab to display the ProDiags Screen.

Click on the **Custom** button to launch the independent ProDiags

This will launch ProDiag Utility.

6 Set idm and idm_send_data tasks as shown below:

1. Click on **Hour** and **Minute** fields of idm_send_data task and select an hour and minute from the pull down lists between 9 AM and 5 PM.

After each selection, click **Yes** on the "Data has been modified. Please confirm it" pop up.

2. Click on **Hour** and **Minute** fields of idm task and select an hour and minute from the pull down lists to schedule 5 minutes before the idm_send_data task.

After each selection, click **Yes** on the "Data has been modified. Please confirm it" pop up.

3. When idm and idm_send_data Hour and Minute selections are complete, click on **Apply** to finalize the schedule updates.

Click **Yes** on the "Are you sure you want to apply changes?" pop up.

Task	Iterations	Day	Hour	Minute	Type	Time Slot
PassChg	1	Daily	05 AM	30	Background	
pd_HouseK...	1	Daily	04 AM	30	Background	
pd_ASC_No...	1	Daily	Hourly	28	Background	
healthpg	1	Daily	04 AM	45	Background	
IIP_HouseK...	1	Daily	06 AM	12	Background	
IIP_HouseK...	1	Daily	12 PM	12	Background	
IIP_HouseK...	1	Daily	06 PM	12	Background	
IIP_HouseK...	1	Daily	12 AM	12	Background	
idm_send_...	1	Daily	07 AM	00	Background	
idm	1	Daily	Hourly	00	Background	

Buttons: New, Delete, Reload, Apply, Help

Status: Reading data... Done

7 There are five new ProDiag tasks that need to be set up according to local hospital system usage. The new tasks are:

- PushTubeUsage
- PushRTS
- sw_download
- hart
- ekg

These tasks need to be added and scheduled (see the following section).

Adding New ProDiag Tasks:

Click New button in the Prodiags utility. This will launch 'Adding New Schedule" pop up. Add the new tasks as indicated in the steps in the figure shown below:

Proactive Diagnostics
Adding New Schedule
Warning: Applet Window

Task	Task Description	Type	Scheduled
idm_send_data	Sends data colle...	Background	Yes
PushRTS	Push RTS Data	Background	Yes
PushRTSAI	Push Complete ...	Background	No
PushScanRx	Push AutoMa D...	Background	No
PushScanUsage	Push Scan Usag...	Background	No
PushTubeUsage	Push Tube Usag...	Background	Yes

Buttons: OK, Cancel

Buttons: New, Delete, Reload, Apply, Help

Status: Reading data... Done

- Click on **New** and the ProDiags utility launches the **Adding New Schedule** pop up.
- Click on the task that needs to be added to the schedule and click on **OK**.
After selection, click **Yes** on the "Are you sure you want to add this task to the schedule?" pop up.
On the "Please edit the newly created scheduled item" pop up, click **OK**.
- Click on **Day**, **Hour** and **Minute** fields of each new task and select a Day (or Daily), an hour and minute from the pull down lists to schedule task.
After each selection, click **Yes** on the "Data has been modified. Please confirm it" pop up.
- When all changes are complete, click on **Apply** to finalize the schedule updates.
Click **Yes** on the "Are you sure you want to apply changes?" pop up.

Recommended default schedules:

PushTubeUsage: Weekly, Wednesday, 10:16 AM
 PushRTS: Daily, 10:07 AM
 sw_download: Weekly, Wednesday, 12:00 PM
 hart: Weekly, Monday, 10:00 AM
 ekg (if available): Weekly, Friday, 4:15 PM

- Close the ProDiags utility browser by selecting **Close** from the upper left browser window pull down menu. The ProDiags browser should close and the InSite™ Interactive Platform Configuration utility should display the HealthPage Tab as shown in Figure 5.

6.3 Configure HealthPage

After completing the ProDiags setup, the IIP Admin Configuration utility automatically switches to the HealthPage.

The Restore State process (after loading Class M software) restored the following information for the HealthPage Reports:

- All e-mail addresses of targeted recipients that were previously loaded
- HealthPage report contents that was previously configured
- HealthPage report schedule that was previously configured

Unless you have recent changes to any of these three items, you can skip the HealthPage Tab of the IIP Admin Configuration utility and select the Device Connection Tab. If you have changes that you want to make (e-mail addresses, report content, report schedule), follow the instructions on the IIP Admin Configuration utility.

The screenshot shows the 'InSite Interactive Platform Configuration' window with the 'Health Page' tab selected. The window has a title bar and a menu bar with 'Information', 'ProDiags', 'Health Page', 'Device Connection', and 'InSite Checkout'. The 'Health Page' section contains instructions for entering e-mail addresses, a 'HealthPage Report Disabled' checkbox, and a text area for 'Enter Address List'. At the bottom, there is a 'Status' bar and buttons for 'DEFAULT', 'CUSTOM', and 'EXIT'. Three numbered callouts are present: 1. Points to the e-mail address input area with the text: 'Enter all GE employee e-mail addresses using the following example: John.Doe@med.ge.com'. 2. Points to the 'HealthPage Report Disabled' checkbox with the text: 'If already not selected, click on the HealthPage Report Disabled box to enable HealthPage reports. Label changes to HealthPage Report Enabled.' 3. Points to the 'Device Connection' tab in the menu bar with the text: 'Click on Device Connection Tab to check connection configuration information.'

HealthPage Screen

1. Enter all GE employee e-mail addresses using the following example:
John.Doe@med.ge.com
2. If already not selected, click on the **HealthPage Report Disabled** box to enable HealthPage reports. Label changes to **HealthPage Report Enabled**.
3. Click on **Device Connection** Tab to check connection configuration information.

6.4 Check Device Connection Setup

When selecting the Device Connection Tab, the IIP Admin Configuration utility automatically displays the Modem Setup Screen.

Make sure that the modem set-up parameters are correct. If not, change the parameters to the correct ones. See Appendix A.

6.5 Check BroadBand Connection Setup

If your system was configured for BroadBand connectivity, the Restore State process (after loading Class M software) restored all the BroadBand setup parameters. If not, enter the proper IP address of the InSite™ Gateway. See Appendix B for instructions.

6.6 Completing the InSite™ Re-Initialization

Follow the instructions provided below to exit the IIP Admin Configuration Utility.

Insite Checkout Screen

1. Click on the Exit button of the IIP Admin Config utility to exit the utility.

2. Click on **OK** on the Exit Confirmation pop-up to exit the IIP Admin Config utility and run the ConfigLink utility.

3. After the shell displays the root prompt, reboot the system to load Insite connectivity parameters. As the system is coming up, you should see the iLinq icon (if the original Checkout included iLinq) as a selectable desktop on the right console monitor.

6.7 Run ConfigLink

Run the InSite™ **ConfigLink** script to populate all configuration information required to fully restore connectivity parameters on the system.


1. Open a shell.
2. Enter root as superuser:
 - a. At the ctuser prompt, type: `su -` <Enter>
 - b. At the Password prompt, type: `#bigguy` <Enter>
3. Change directory to InSite™. At the root prompt, type: `cd /usr/g/insite` <Enter>
4. Run the ConfigLink script. At the root prompt, type: `ConfigLink` <Enter>

Wait for the shell to return the root prompt to ensure the **ConfigLink** process completes. When complete, exit root and close the shell.

6.8 Turn OFF telnet in PNF

The default setting in PNF allows telnet from all IP addresses. To make the system more secure, turn off telnet via the commands below. If telnet is required for systems within the hospital network, individual IP addresses can be added as needed.

1. Open Common Service Desktop
2. Click on Configuration Tab
3. Click on Configure PNF.
4. Follow steps 1-4 shown in picture below.

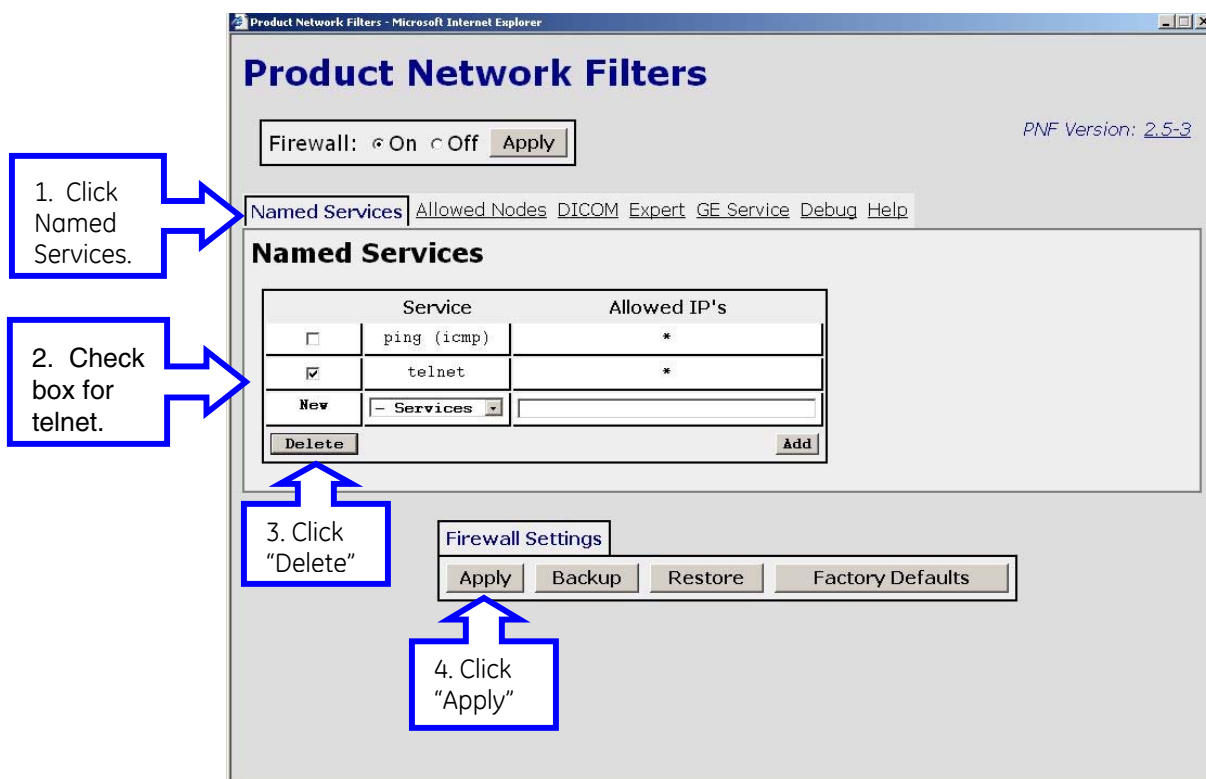

6.9 Reboot System

1. Click on the pink Shutdown button.
2. From the Shutdown selection pop-up, click on RESTART, then OK, to reboot the system.

6.10 Save System State

At this time, save Reconfig portion of the system state using the Save System State utility.

7 InSite™ Update Check List/Feature Demo to Customer

7.1 Check List

Field Engineer Actions

The field engineer is to demonstrate the InSite™ Update (Software Download/Update Management) feature to the radiology center's manager, or supervisor, to determine how their system is to be configured:

After the LFC is completed:

- Demonstrate InSite™ Update (download and install) to the InSite™ broadband connected customer (see 7.2).
–or–
- Walk-through the InSite™ Update marketing material (kit p/n 5164884) with the customer to ensure understanding and value InSite™ Update delivers.

Important Discussion Points

The field engineer must ensure the customer understands that the InSite™ Update feature delivers the following:

- Access to the latest system updates as soon as they are available for their specific hardware and software configuration.
- No need to schedule an FE to receive the update.
- The customer determines the timing for update installations which deliver:
 - o Increased system uptime
 - o Improved system performance/quality

Options:

After completing the InSite™ Update demo, (or discussion of the marketing material), if the customer is concerned about the feature, the feature can be temporarily disabled and followed up with a live demo when next update is available.

When next time when an update is available, schedule time with the customer and demonstrate the feature by manually down loading the update and then auto installation. If the customer does not want the feature enabled at that time, then explain to the customer that there are two aspects of the feature, Auto-Download and Auto-Install. Auto-Download can be enabled without the need for Auto_Install thereby requiring FE to install the update.

General Notes:

- The system must be broadband connected to maximize the full benefit of the automatic InSite™ Update feature.

7.2 Demo 'InSite™ Update Feature to the Customer

Installation of demo software and display to customer


- Open Unix shell
 - Type: `su -` <Enter>
 - Type password: `#bigguy` <Enter>
 - Type: `swd_demo -e` <Enter>

(This command sets up an update demo that enables auto SW downloads)

The following message will appear:

Successfully set Demo mode of SW download!

- Reboot System
- On Boot-Up, a pop-up menu will ask the user to Install the Demo update Now or Later – see screen shot below

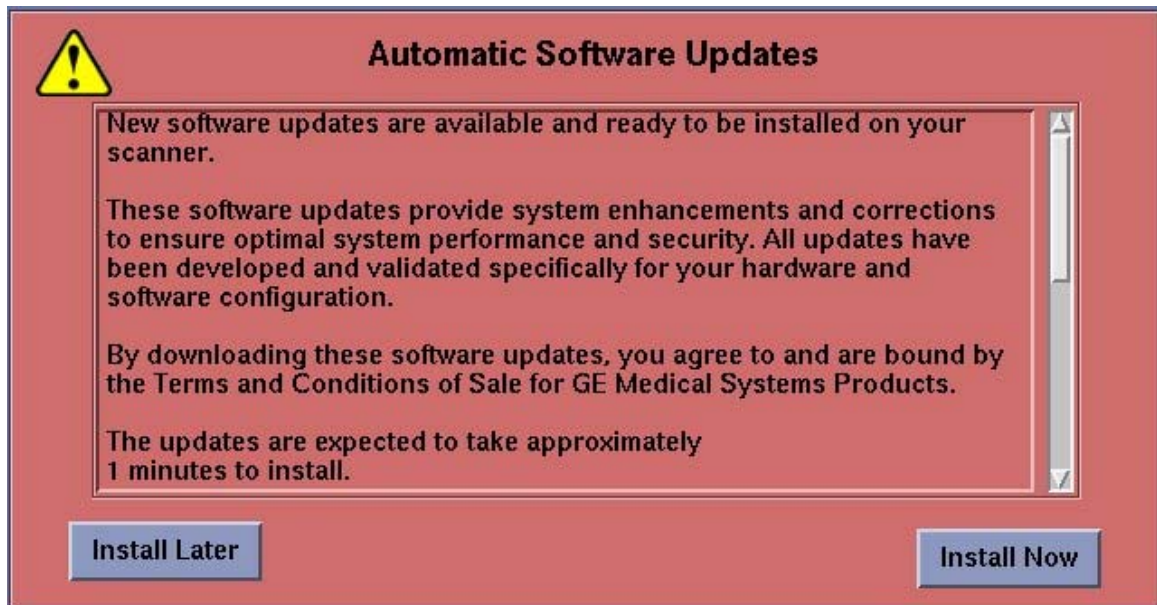

- Select Install Now (If Install Later is selected, this window will pop up at each reboot)
The Demo software update will be installed in the scanner and the following message will be displayed



- To verify that the update has been installed: Go to common serviced desktop and select *Utilities/Update Management/Update Status* to view the following screen:

Update_Status			
Update Name	Status	Action	Install Time
SWD-DEMO	Installed on the scanner	NONE	

Clean up the demo software

Clean up demo software after the demo to the customer:

- Open Unix shell
- Type: `su -` <Enter>
- Type password: `#bigguy` <Enter>
- Type: `swd_demo -c` <Enter>

The following message will appear:

Successfully clean up Demo environment of SW download!

Do you need to disable auto install on this scanner? [y/n]

Answer "n" (no), unless the customer indicates that they do not want this feature enabled.

8 Finishing Up

Leave the following documents with the customer:

- Advanced Applications on OC OP Info, 2404359-247
- CT Light Speed series CD, 2369740-247
- New Features Guide for 04MW13.7 and 04MW44.9 Software, 5166132-100
- Application Tips and Precautions, 5167065-100
- “Customer Letter,” 5142333-147 (the correct language version).
- Marketing envelope, 5164884 which contains a tri-fold pamphlet, Tent Communication and an InSite™ Update mouse pad. This applies only to InSite™ Entitlement systems.

Remove the items from the envelope. Place the Tent Communication card on the keyboard. Leave the mouse pad and the leaflet on the console.



9 FMI Closure


Close FMI per regional procedure and close.

Remove FMI Instructions from the customer site.

APPENDIX A

Modem Setup

When selecting the Device Connection Tab, the IIP Admin Configuration utility automatically displays the Modem Setup Screen as shown below.

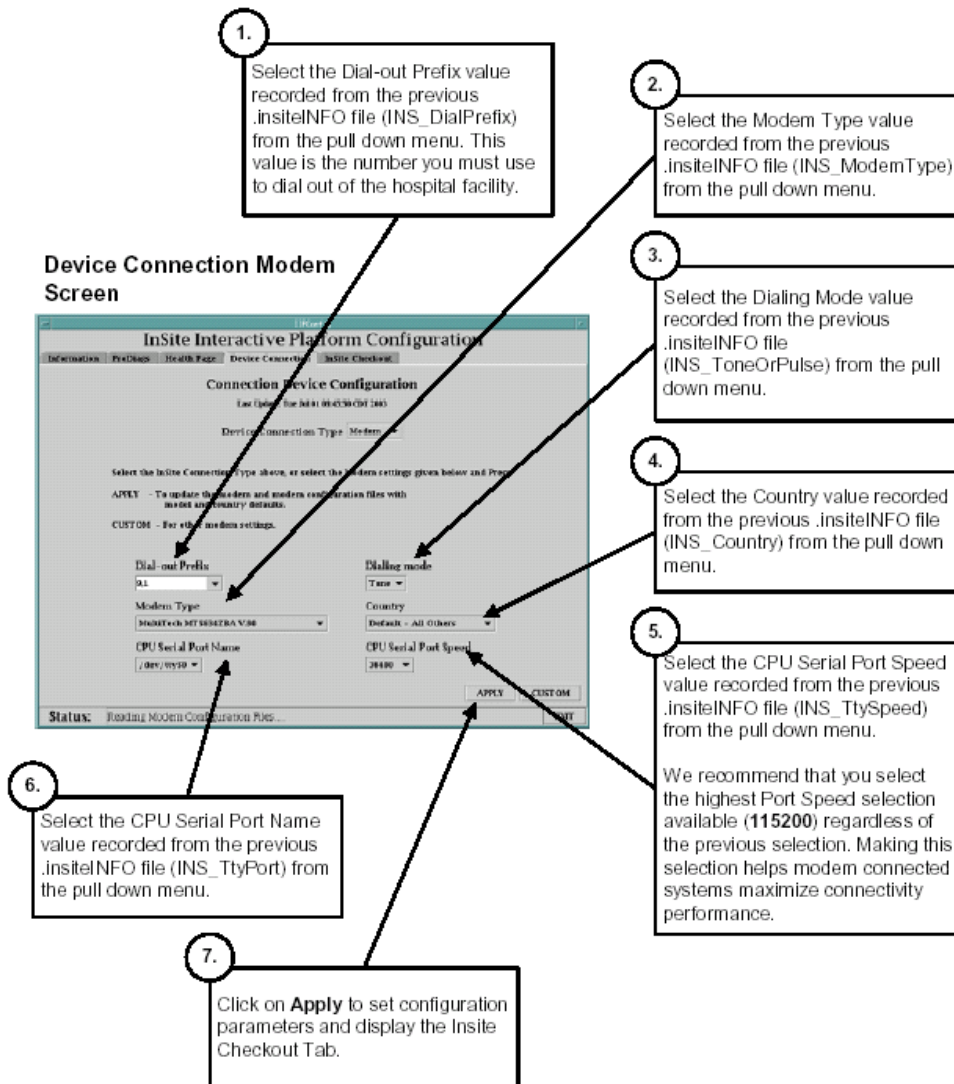

If your system was configured for Modem connectivity, the Restore State process (after loading Class M software) restored all the Modem setup parameters. Verify with the InSite™ information that you recorded before starting the loading process.

- If there were no Modem parameters in the .InSite™INFO file but there were BroadBand parameters, skip to the Check BroadBand Device Connection Setup procedure.
- If there were Modem parameters in the .InSite™INFO file but the parameters displayed in the Modem Setup Screen do not match what you previously recorded, enter the correct values as required.

NOTE: The Modem configuration parameters should be pre-loaded from the Restore State process. Do not change settings unless these values were not properly restored.

Appendix B

Broadband Connection Set Up

If your system was previously BroadBand InSite™ connected, click on the Network button of the IIP Admin Configuration utility Device Connection Tab to display the Network Setup Screen.

If your system was configured for BroadBand connectivity, the Restore State process (after loading Class M software) restored all the BroadBand setup parameters. If not, enter the proper IP address of the InSite™ Gateway as shown in the figure below.

NOTE: The BroadBand Network configuration parameters should be pre-loaded from the Restore State process. However, the on-site FE who performed the original InSite™ Checkout of the system may have used the hospital's imaging LAN default gateway IP Address (selection was available in the previous release of IIP).

DO NOT USE THE HOSPITAL WORKING LAN DEFAULT GATEWAY IP ADDRESS AS THE InSite™ GATEWAY IP ADDRESS OTHER THAN AS A TEMPORARY MEANS TO CONNECT TO THE ON LINE CENTER!! Contact the hospital network administrator to determine the right gateway address.

Device Connection Network Screen

1. Enter the IP Address of the Insite Gateway to GE Online Center as recorded from the .insiteINFO file (INS_GatewayIP). Do not enter leading "zero's" into the four fields.
2. Click on **Apply** to have the IIP Admin Config utility ping the entered IP Address, load the network parameter, and display the Insite Checkout Tab.
3. Click on **OK** on Checkout pop-up

FMI	FIELD MODIFICATION INSTRUCTION	FMI
25365	06MW03.5 SW Update for CT-LS 1X/2X/3X/4X/5X & HS QX/i with GRE Consoles	PAGE 24 OF 24

This Page Intentionally Blank.