

**ADMINISTRATIVE INFORMATION****FMI No: 25367****FMI Title: Helios 2214768 & 2214768 -2 Collimator Xray Scatter Filler Bar for PET/CT & CT Collimators**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS - AM:	Issue Date:
☒ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked: See "Effectivity" Section****FMI Instructions Part No.: 5142977-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25367

FMI Title: Helios 2214768 & -2 Collimator Xray Scatter Filler Bar

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments _____

Effectivity

The following systems are part of FMI 25367's effectivity:

All CT and PET/CT systems equipped with the Helios Collimator, tracked model numbers:

- 2214768 Helios Collimator and
- 2214768-2 Helios 2 Collimator.

Purpose

This FMI is to install a small filler bar (p/n 5134695) on all Helios collimators, 2214768 and 2214768-2.

The filler bar is to be installed on the CT collimator frame to reduce the risk of potential x-ray containment issues at the seam between the filter assembly and the collimator frame of the 2214768 and 2214768-2 collimators.

Prerequisite FMIs


None

Related FMIs

This FMI is not batched with any other FMI.

This FMI is related / similar to FMI 25360 that deployed the collimator filler bar to all Helios –3 and –4 gantries.

Furnished Materials

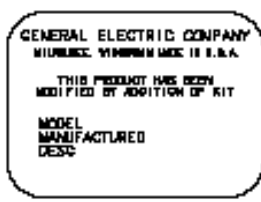FMI 25367, Kit P/N 5143925			
Qty	Description	Part #	Comments
1	FMI Document	5142977-100	
1	M4 x 0.7 bottoming-style hand tap	5135684	
1	Rating Plate - Sticker	46-302200P12	
1	Notched Filler Bar Assembly	5134695	
1	3mm Hex Shaft Tool Bit	5151374	McMaster Carr p/n 7397A43 Fits torque driver, p/n 46-282782P1

Tool Required

2-36 in-lb torque driver (p/n 46-282782P1) from service tool kit (p/n 46-268445G1)

Personnel Requirements

1 FE

1 hour, labor

1 hour, travel

Procedure

Refer to the Service Methods CDROM for detailed instructions in support of these FMI tasks (e.g., LOTO).

This procedure provides instruction to install the collimator filler bar to 2214768 and 2214768-2 collimators.

1. Remove gantry covers to gain access to the CT collimator.
2. Visually inspect the CT collimator assembly to determine if the collimator filler bar has already been installed.

Refer to Figure 3 & 4.

3. If the collimator filler bar is installed, STOP. Do not replace it. Close the FMI with completion code = 1, "Modification Complete or previously modified."

If the collimator filler bar is not installed, continue executing this FMI.

4. Turn off the gantry axial enable and HVDC enable switches.
5. Turn off the gantry 120 Vac switches to completely power-down the gantry.
6. Perform Lockout Tagout on the gantry (LOTO).
7. Rotate the gantry so that the CT collimator is at the 6 o'clock position.
8. As a precautionary measure, cover the CT collimator output port with a towel to prevent debris from getting onto the port.

9. Using the 3mm tool bit (included, p/n 5151374) and torque driver (p/n 46-282782P1) from the service tool kit (p/n 46-268445G1), remove the two (2) existing M4 x 20mm screws from the front of the collimator.

Refer to Figure 1 & Figure 2 for the location of these screws; discard them.



Figure 1 - Example of screw being removed

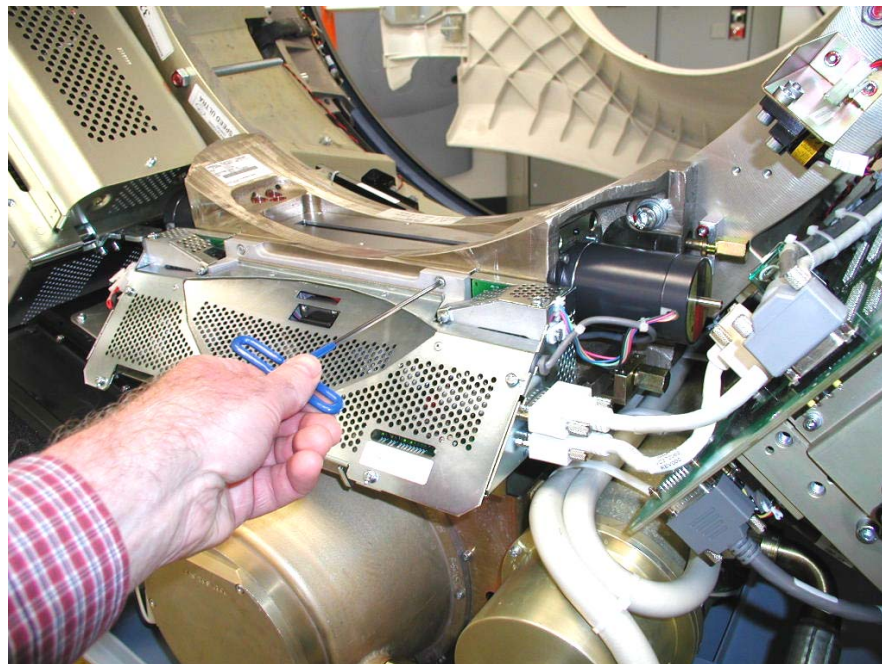


Figure 2 - Example of screw being removed

10. Using the hand tap (included, p/n 5135684), clean these screw holes of LocTite by gently turning the tap in each screw hole until it lightly bottoms out.

To avoid cross-threading, or breaking the tap, do not force.

11. With the notch of the filler bar (included, p/n 5134695) facing downwards, attach the filler bar to the collimator assembly.

Refer to Figure 3 & Figure 4.

NOTE: The filler bar only fits in one direction – notch downward as shown in Figure 4.

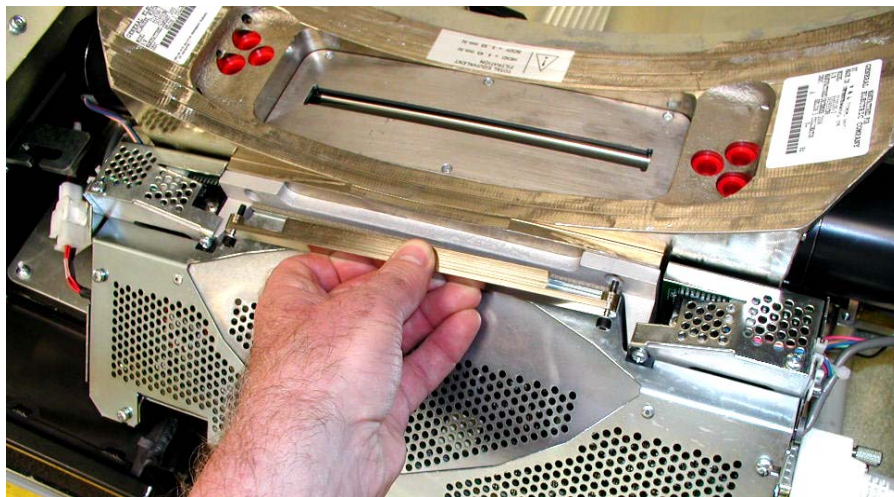


Figure 3 - Installing filler bar

12. Tighten the two M4 screws to 2.3 N-m, or 20.4 lb-in.

Use the 2-36 in-lb torque driver (p/n 46-282782P1) and the 3mm tool bit (included, p/n 5151374) to properly torque the screws.

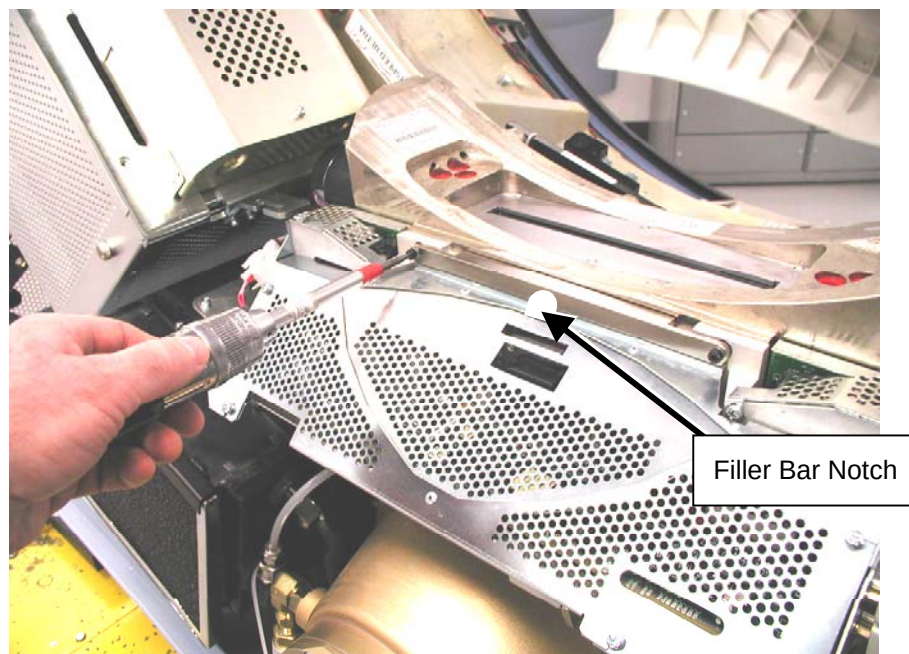


Figure 4 - Torque the M4 Screw Using Torque Driver, 46-282782P1

13. Place the rating plate sticker (included, 46-302200P12) on the collimator frame as shown in the picture, Figure 5.

Do not cover up any text on the existing labels.



Figure 5 - Rating Plate Sticker Placement

Finishing Up

14. Remove LOTO per detailed instructions in the Service Methods CDROM.
15. Turn on 120 Vac power to gantry.
16. Turn on the gantry axial enable and HVDC enable switches.
17. Replace gantry covers.
18. With covers in place, rotate the gantry to ensure there is no mechanical interference between the collimator filler bar and gantry cover components. Adjust covers, if required (no cover adjustment expected).
19. Perform check scan (for example, FastCal), to confirm system performance before returning system for patient scanning.
20. Close FMI service dispatch per local procedures.