

**ADMINISTRATIVE INFORMATION****FMI No: 25378****FMI Title: HP60 Gantry Aux Box Mounting Bracket Torque Check for CT & PET/CT HP60 Gantries**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support		
GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked: See "Effectivity" Section****FMI Instructions Part No.: 5162550-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25378

FMI Title: HP60 Gantry Aux Box Mounting Bracket Torque Check
for CT & PET/CT HP60 Systems

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments

Effectivity

CT & PET/CT systems based on the HP60 gantry with tracked model numbers listed below and manufactured after 15-Aug-05 (specific effectivity was generated based on CT Manufacturing—Waukesha data):

- **2383550-11**, HP60 LIGHTSPEED5.0J GDAS4 EMC ED2 ASTEC
- **2383550-70**, HP60 GDAS4 EMC ED2 ASTEC FOR PET
- **2383549-11**, HP60 LIGHTSPEED5.0J GDAS8 EMC ED2 ASTEC
- **2383549-70**, HP60 GDAS8 EMC ED2 ASTEC FOR CT (for PET)
- **2378661-11**, HP60 LIGHTSPEED5.0J GDAS16 ED2 ASTEC

Purpose

This FMI is to confirm the torque on five (qty 5) bolts that are part of the HP60 Aux Box mounting bracket, p/n 2370279.

Personnel Requirements

1 FE Resource
 1.0 hour Labor
 1.0 hour Travel

Furnished Materials

FMI 25378 Kit P/N = 5162552			
Qty	Part #	Description	Comments
1	5162550-100	FMI Instruction	
1	5155306	Torque Bit, 8mm hex, ball, 3/8" socket drive	Machinery & Factory p/n FABM8E

Tools & Materials Required

Qty	Part #	Description	Comments
1	46-315668P2	Torque Wrench, 3/8" drive, 200-1200 in-lb	Found in FE Service Tool Kit, p/n 46-268445G1
1	46-170686P3	Loctite 242	Found in HP60 Installation Support Kit, p/n 2319891-3
1	2284015	Torque Bit, 10mm hex, 3/8" socket drive	Found in HP60 Installation Support Kit, p/n 2319891-3
1	2284016	Torque Bit, 5mm hex, 3/8" socket drive	Found in HP60 Installation Support Kit, p/n 2319891-3
1	N/A	3/8" drive extender, 6"	

Procedure

This instruction details the process to confirm the torque on five (5) M10 bolts that are part of the CT and PET/CT HP60 Aux Box mounting bracket, p/n 2370279. Refer to Figure 1: Five Bolts on Aux Box Mounting Bracket to Torque, for an illustration of the five bolts whose torque must be checked.

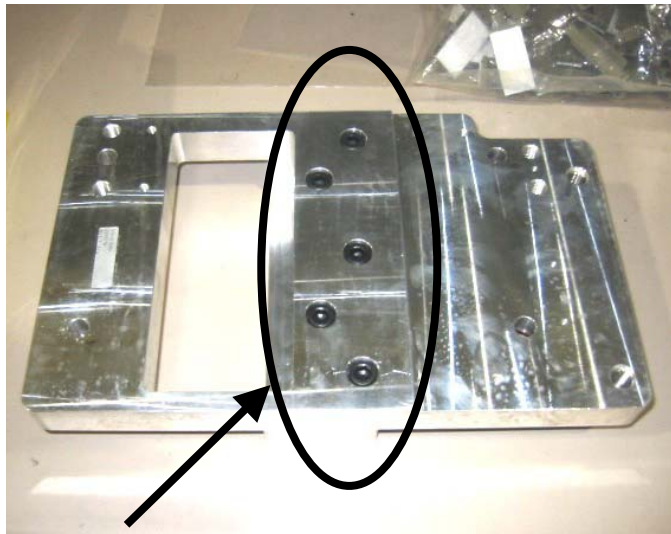


Figure 1: Five Bolts on Aux Box Mounting Bracket to Torque

The Aux Box mounting bracket secures the Aux Box to the rotating portion of the HP60 gantry as illustrated in Figure 3: HP 60 Aux Box Subassembly. The following steps detail the part-removal process to gain access to the Aux Box mounting bracket.

- 1) Remove gantry covers to gain access to the Aux Box and gantry service switch control board.
 - a. right-side gantry cover,
 - b. left-side gantry cover, and
 - c. left-top gantry cover.
- 2) Turn off the power to the gantry and perform lock-out / tag-out.
- 3) From the left side of the gantry, rotate the gantry by hand to position the Aux Box at the 10 o'clock position to make it most accessible.
- 4) Turn the gantry lock to engage the ring, preventing it from rotating.

Refer to Figure 2: Gantry Lock Mechanism.

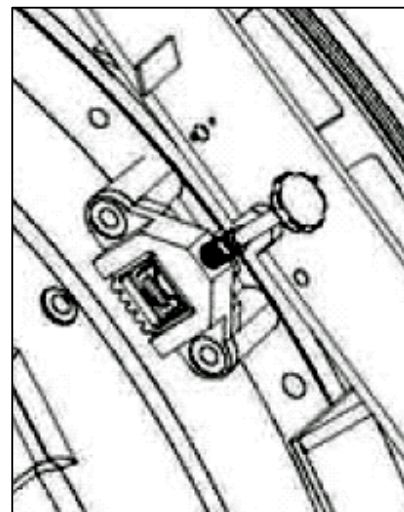


Figure 2: Gantry Lock Mechanism

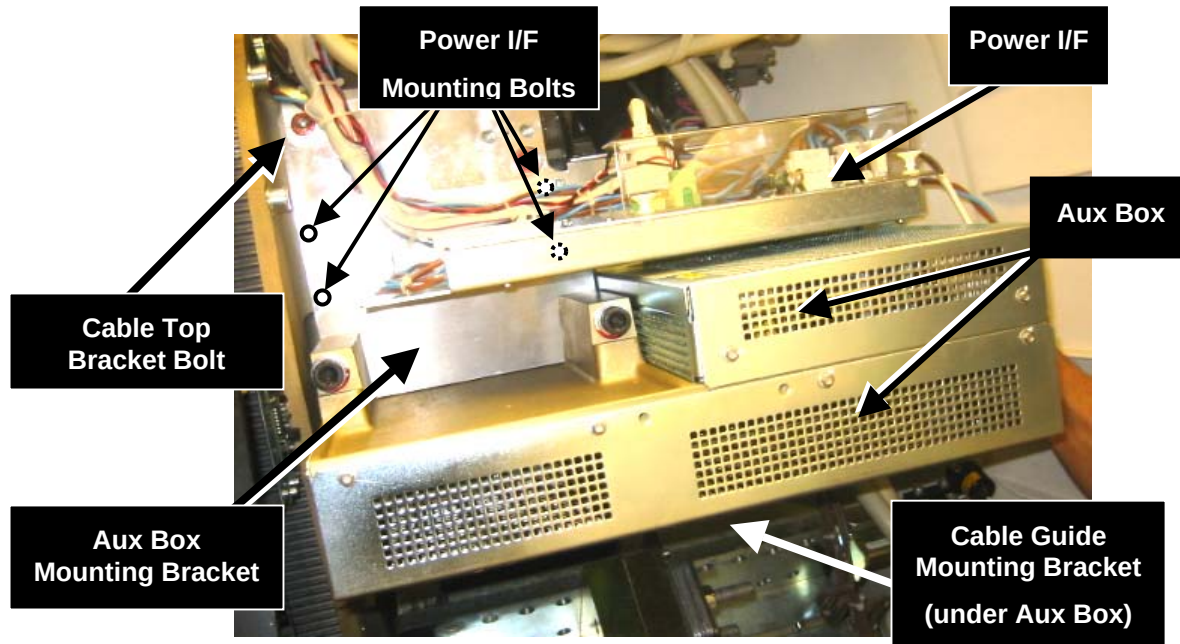


Figure 3: HP 60 Aux Box Subassembly

- 5) Remove the two (2) M12 x 150 Aux Box mounting bolts – use 10mm bit from the HP60 Installation Support Kit, p/n 2319891-3.

Refer to Figure 3: HP 60 Aux Box Subassembly.

- 6) Move the Aux Box to the side (rest it on the gantry's front cover bracket) to gain access to the Aux Box mounting bracket.

NOTE: There should be no need to remove any cable connections to the Aux Box, but if your cables are too tight, remove cables to allow the Aux Box to be removed.

- 7) Remove one (1) M12 x 30 bolt which mounts the cable top bracket (p/n 5111277) to the Aux Box mounting bracket – use 10mm bit.

Refer to Figure 3: HP 60 Aux Box Subassembly.

- 8) Remove four (4) M6 x 12 Power I/F mounting bolts – use 5mm bit from the HP60 Installation Support Kit, p/n 2319891-3.

Refer to Figure 3: HP 60 Aux Box Subassembly.

- 9) Remove two (2) M6 x 12 cable guide mounting bracket (p/n 2370277) bolts – use 5mm bit.

Refer to Figure 3: HP 60 Aux Box Subassembly.

- 10) Remove four (4) M12 x 40 mounting bolts holding the Aux Box mounting bracket to the rotating gantry – use 10mm bit.

Refer to Figure 4: Aux Box Mounting Bracket Removal

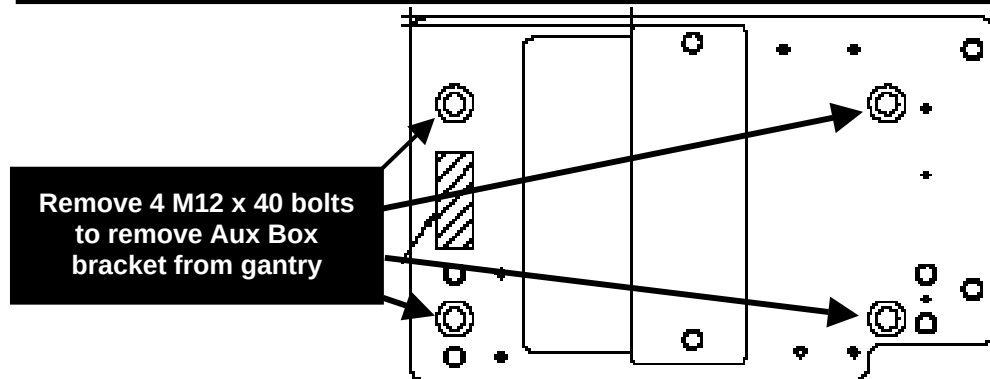


Figure 4: Aux Box Mounting Bracket Removal (4 Bolts)

- 11) Remove the Aux Box mounting bracket from the gantry.
- 12) Retrieve the torque wrench (p/n, 46-315668P2) from your service / torque tool kit (p/n 46-268445G1) and the 8mm hex torque bit provided in this kit.
- 13) Place the Aux Box mounting bracket flat on the floor and torque the five (5) M10 bolts to 25 ft-lb (33.8 N-m).

IMPORTANT: In Figure 5, capture the degrees swept by the torque wrench to tighten each bolt.

Report findings (FE Name, SID #, degrees swept to tighten each bolt) by calling:

- Gerry Bass (OLCE) @ 262-896-2867 or
- Steve Surprenant / Randy Ulery (OLCE) @ 262-544-3081

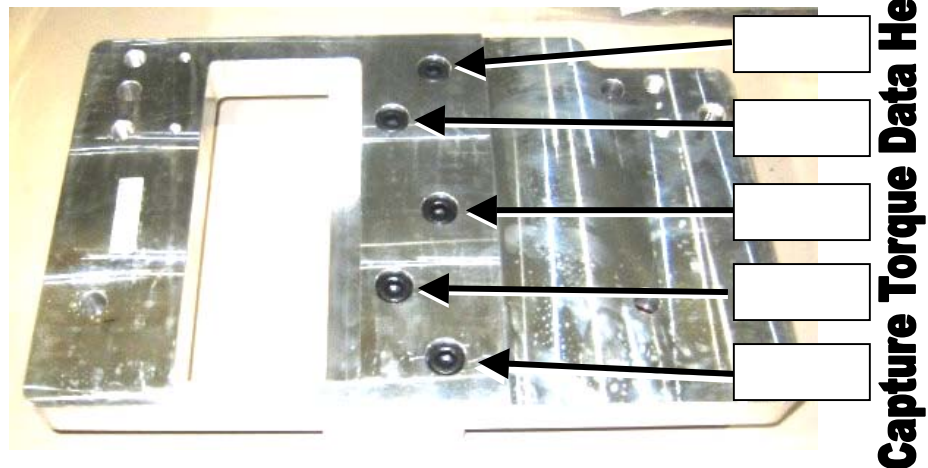


Figure 5: Degrees Swept to Torque Each Bolt

- 14) Retrieve the bottle of LocTite 242 (p/n 46-170686P3) from the site's HP60 Installation Support Kit (p/n [2319891-3](#)).
- 15) Re-mount the Aux Box mounting bracket to the rotating gantry; apply LocTite 242 to the four (4) M12 x 40 mounting bolts and torque to 49 ft-lbs (66 N-m).

NOTE: Apply torque guidelines; tighten all bolts to ½ of the torque specified, then complete to full-torque specified.

- 16) Apply Loctite 242 to one (1) M12 x 30 cable top bracket mounting bolt and torque to 49 ft lbs (66 N-m).
- 17) Re-mount the Aux Box and torque two (2) M12 x 150 bolts to 49 ft-lbs (66 N-m).
- 18) Re-mount the Power I/F assembly; apply Loctite 242 to the four (4) M6 x 12 mounting bolts and torque to 5.9 ft-lbs (8 N-m).
- 19) Re-mount cable guide bracket; apply Loctite 242 to the two (2) M6 x 12 mounting bolts and torque to 5.9 ft lbs (8 N-m).
- 20) Check for any loose cable connections to the Aux Box (if removed, reconnect cables) and replace any removed ty-wraps. Refer to Figure 6: Aux Box Cable Connections.
- 21) Release the gantry lock mechanism.
- 22) Rotate gantry by hand to make sure no components are loose or interferences exist.
- 23) Restore power and rotate gantry 2-4 revolutions using service switch.

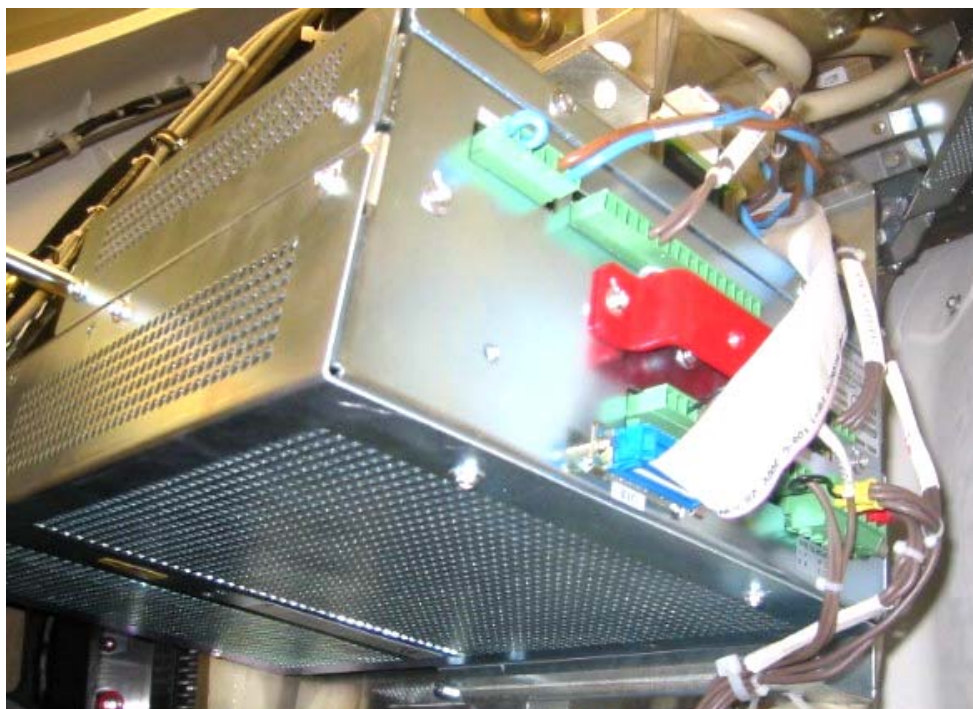


Figure 6: Aux Box Cable Connections

Finishing Up

- 24) Replace gantry covers.
- 25) Perform check scan (for example, FastCal) to confirm system performance before returning system for patient scanning.
- 26) Close your FMI dispatch.
- 27) Confirm: Call OLCE representatives listed in Step #13 and leave your torque results:
 - Gerry Bass (OLCE) @ 262-896-2867 – or –
 - Steve Surprenant / Randy Ulery (OLCE) @ 262-544-3081.