

**ADMINISTRATIVE INFORMATION****FMI No: 25389****FMI Title:** 07MW11.10 Software Quality Release for LightSpeed 1X, 2X, 3X, 4X, 5X (not RT16 or Xtra), HiSpeed QX/i, and BrightSpeed 4, 8, and 16 slice (not BrightSpeed Select).

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS – AM:	Issue Date:
☐ Mandatory	GEMS – E:	
☐ Mandatory on Request	GEMS – A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS – Americas	GEMS - Europe	GEMS – Asia
GE Medical Systems On-Line Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**PAY SPECIAL ATTENTION to the**

- “Prerequisites” and
- “Section 3.1, Need to Know Before Performing LFC”

sections **BEFORE performing the load-from-cold** of software release 07MW11.10 provided in this kit.**Unit Tracked:** Multiple Host Console Model Numbers**FMI Instructions Part No.:** 5191771-100, Rev 4**Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS – Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25389

FMI Title: 07MW11.10 Software Quality Release for LightSpeed 1X, 2X, 3X, 4X, 5X (not RT16 or Xtra), HiSpeed QX/i, and BrightSpeed 4, 8, and 16 slice (not BrightSpeed Select)

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments _____

Effectivity

This FMI is to be installed on CT-only LightSpeed 1X, 2X, 3X, 4X, 5X, RT4 (not LightSpeed RT16 and not LightSpeed Xtra), BrightSpeed 4, 8 and 16 slice (not BrightSpeed Select), and HiSpeed QX/i systems.

The CT system must be configured with a GRE (GOC3 or GOC4) host computer console with one of the following tracked model numbers:

• 2377708	GOC3
• 2377708-4	
• 2377708-7	
• 2377708-10	
• 2377708-13	
• 2377708-20	
• 2377708-21	GOC3
• 2377708-22	GOC4
• 2377708-30	
• 2377708-31	
• 2377708-32	
• 2377708-40	
• 2377708-41	
• 2377708-42	
• 2377708-43	
• 2377708-50	
• 2377708-51	
• 2377708-52	
• 2377708-53	
• 2377708-120	
• 5111750	GOC4
• 2377708-520	GOC4 (Jarrell)
• 2377708-530	
• 2377708-540	
• 2377708-550	
• 2377708-551	
• 2377708-552	
• 2377708-553	
• 5377708-120	
• 5377708-520	GOC4 (Jarrell)

DO NOT INSTALL this FMI on any:

- PET/CT,
- LightSpeed RT16,
- LightSpeed Xtra, or
- BrightSpeed Select

(BrightSpeed Select has MX165 glass tube, Insert Model # 46-309300G3, Casing Model # 2199077-7, & does not have a HEMIT tank).

Purpose

This FMI delivers software release 07MW11.10 to CT-only systems to improve the quality of LightSpeed 1x2x3x4x5x, BrightSpeed 4/8/16 slice, and HiSpeed QX/i systems. No new features are delivered in this software release.

Software release 07MW11.10 also delivers software to support the transition to a different Hewlett Packard (HP) PC graphics card due to end of life (FX1400 transitioning to FX1500).

For systems with PS/2-cabled trackballs, this FMI also delivers a PS/2-to-USB cable converter to ensure the trackball functions with the operating system, Linux OS Rls 6.2.9, installed with this software release 07MW11.10.

Prerequisites

The CT-only system must be loaded with one of the software releases listed below before installing this FMI's software:

- 06MW03.5, FMI 25365, kit part number 5137245 -- or --
- 06MW29.7, software collector part number 5193989.

Check the software release currently loaded by:

1. Entering Service/Common Service Desktop.
Double-click "Service" icon on right monitor of host console.
2. View "Home" tab, "Current System Status" area.
3. View "Application Software" field, "Information" column.

Look for "06MW03.5" or "06MW29.7" at the beginning of the string (or "07MW11.10;" see "**STOP**" note below).

STOP if your system does not have one of these software releases already loaded. One of the software releases listed above must first be installed before proceeding with this FMI.

STOP if your system is a PET/CT system. This FMI is not to be installed. Close this FMI dispatch as Code 2 and return the FMI kit.

07MW11.10 SOFTWARE ALREADY LOADED?

STOP -- If your system is already running the 07MW11.10 software release, do not execute this FMI's load-from-cold. Jump to Section 8 entitled "Service Pack 1.5 Software Load" and load Service Pack 1.5 software.

Furnished Materials

The following items are included in FMI 25389, kit part number 5191771, Rev 3:

Item	Name	Rev	Qty	Description
1	5191771-100	3	1	FMI 25389 Instruction - 07MW11.10
2	5167277-200	6	1	Service Information CD (LFC), FX_Xtream SW & Upgrade Procedures
3	5220443	3	1	07MW11.x Software Collector with Rating Plate for HS QX/i, LS 1x2x3x4x, LS 5x 4,8,16slice, and RT4
4	5234959-15	1	1	Service Pack 1.5 CDROM for sw rls 07MW11.10, 07MW18.4, 07BW8.3, and 07BW8.4 w/ Product Locator
5	5191796-3	1	1	BIOS 2.10 Floppy Set for HP 8200 only, BIOS Floppy and Setup Floppy
6	5138295	1	1	BIOS 1.13 Floppy Disk for HP 8000 only
7	5195935	1	1	Trackball USB-to-PS2 converter. Fixed for CTT 6.x Operating System.
8	2337414	2	1	DVD RAM MEDIA, 9.4 GB
9	2369740-247	9	1	CT LIGHTSPEED SERIES CD
10	5119787-100	8	1	LS Series QuickStep Guide -- English
11	5119787-247	8	1	LS QuickStep Guide
12	5141172-247	3	1	BrightSpeed LRG Operator Manual
13	5159870-100	3	1	BrightSpeed QuickStep_English
14	5159870-247	3	1	BrightSpeed QuickStep Guide - Multi-language CD
15	5198350-100	1	1	2007 Common Platform Software Release Application Tips and Workarounds, Paper Version
16	5198350-247	1	1	2007 Common Platform Software Release Application Tips and Workarounds, Multi-Language CD
17	5231794-148	1	1	Product Matrix Manual
18	5119787-141	8	1	LS Series Quick Step Guide -- China
19	5159870-140	3	1	BrightSpeed QuickStep_Japanese
20	5159870-141	3	1	BrightSpeed QuickStep_Chinese
21	5166779-141	2	1	BrightSpeed SFDA Document
22	5183783-247	1	1	CT--PET Advanced Application II on O OP Info-multi-language collection
23	5198350-140	1	1	2007 Common Platform Software Release Application Tips and Workarounds, Paper Version-Japanese
24	5198350-141	1	1	2007 Common Platform Software Release Application Tips and Workarounds, Paper Version-Chinese
25	5252415-199	1	1	Learning and Reference Guide Addendum- English, French, German, Italian, Portuguese, Spanish, Japanese, Chinese and Korean language
26	5119787-140	8	1	LS Series Quick Step Guide -- Japan

Tools Required

- 07MW11.10 Class C Software, p/n 5220446
- 07MW11.10 Class M Software, p/n 5220446-2

NOTE: 07MW11.10 Class C and M CDRoms were distributed to FEs via a TAB #182 distribution in advance of this FMI. If you have not received them, contact your process leader to ensure you are on TAB #182's distribution list. Order them through the GEMS CoakleyTech website at <http://coakleytech.med.ge.com/cgi-win/default.aspx>

Personnel Requirements

Travel:	1 Hour
Labor:	6 Hours
Field Engineer:	1

Product Safety Fixes

- System shutdown due to XIO error (multiple PSR's)
- Scan aborts, start scan fails to light, table or cradle not ready due to firmware start up issues on LightSpeed 5.X and BrightSpeed systems (multiple PSR's)
- Table moving down on it's own
- Advanced Vessel Analysis gives incorrect measurements if adjust tracking
- Mixed languages found in VoxTool / Reformat (fixed by Service Pack 1.5)

Product Quality Fixes


- mA limited if Tube Warm Up skipped on LightSpeed Pro16 systems (Hercules tube)
- Warning message posted that Auto mA is turned "off" when patient orientation and patient entry do not match the scout series
- SmartPrep reliability improvements
- SmartView reliability improvements
- Start up Shut down reliability improvements
- Unable to start New Patient (pixMgmt crash)
- Images fail to Auto View on LightSpeed 4.X 3 IG systems if SmartPrep used
- Edit Patient fails due to "EditPatientData.ct" file not updated correctly
- Alt-F4 crashes Display process or Display Manager process

- Time drift and time jump
- Daylight Savings Time, now beginning 2nd Sunday in March and ending 1st Sunday in November) has been implemented
- Multi language Auto voice defaults to reconfig default voice at End Exam
- Images fail to transfer to a Kodak PACS
- Images fail to copy or restore from CD/DVD
- Unable to Save Multi Step protocols in Reformat or Volume Viewer

Enhancements / Improvements

- Auto Transfer by Series for Scouts
- Auto Transfer for Dose Report
- Prospective Multi Recon for Cardiac scan type
- User selectable series description for Retro recon, Reformat and Volume Viewer
- IQ Snap -- reserve the scan data files related to images with IQ issues for further investigation
- Full and Partial Anonymous patient levels
- Scan date added to Large Image Font
- Screen Save of Exam or Series text page
- Patient Name and Patient ID added to update parameters in Connect Pro
- Save / Restore System State time optimized (time reduced)

Process Flow

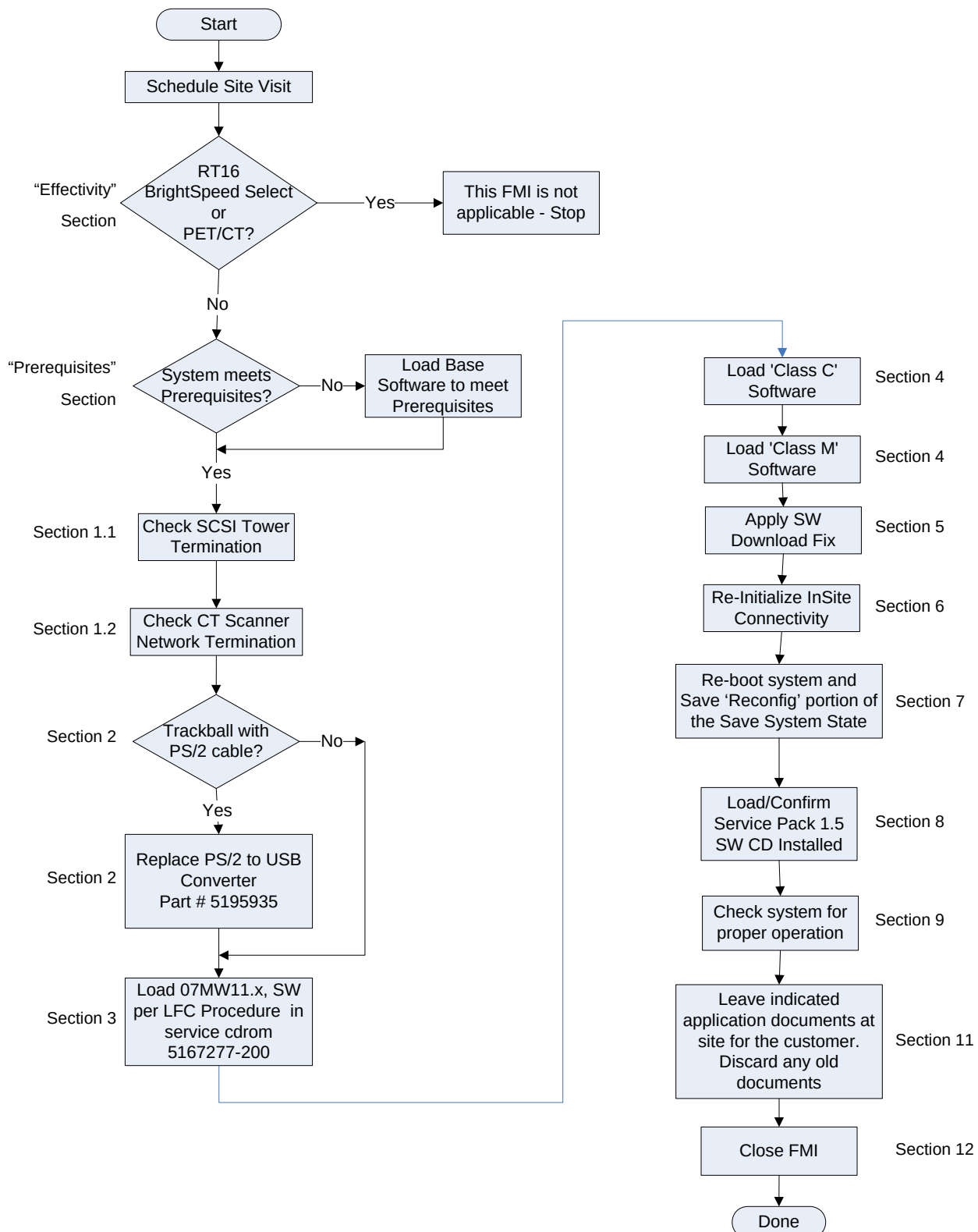

Procedure

1. Preliminary Checks

1.1 SCSI Tower Check

Check your SCSI tower to ensure its SCSI bus is not doubly terminated.

Certain SCSI Towers had a problem with double termination that affected SCSI function, such as Save System State. It is important to check your SCSI Tower to see if yours is affected.

1. Check the part/model number on the bottom of your SCSI tower.
 - a) If the SCSI tower part number is p/n 2351055 or p/n 5167410, then skip this section and proceed to Section 1.2 entitled "CT Scanner Network Check."
 - b) If the SCSI tower p/n is 2368838, then execute Service Note 06080312a.

This service note is included in the Appendix of this FMI instruction.

The service note may be accessed at the following link:

http://www.ct.med.ge.com/tech_pubs/2360027/sn/SN06080312a.pdf

1.2 CT Scanner Network Check

Check the CT scanner's network (host console to gantry & table) / media cable converter to ensure it is correctly terminated to 50 Ohms.

The network media converter, p/n 2377227 (Milan model MIL-C1111), has been suspected of causing communication issues that may cause a scan not to complete correctly.

The following network media converter setup information is specific to the Milan model MIL-C1111, p/n 2377227.

2. Locate the Network Media Converter; refer to **Illustration 1**.
3. Set / check its DIP-switch settings per **Illustration 2** and ensure no BNC terminator is attached (enable 50 Ohm termination at end of network).
4. Set / check its MDI-X/MDX slide switch setting per **Illustration 3** (enable connecting to workstation / another network converter).

Illustration 1: Network Media Converter Location, Affixed to Side of Chassis (Network Media Converter p/n 2377227, Milan model MIL-C1111)



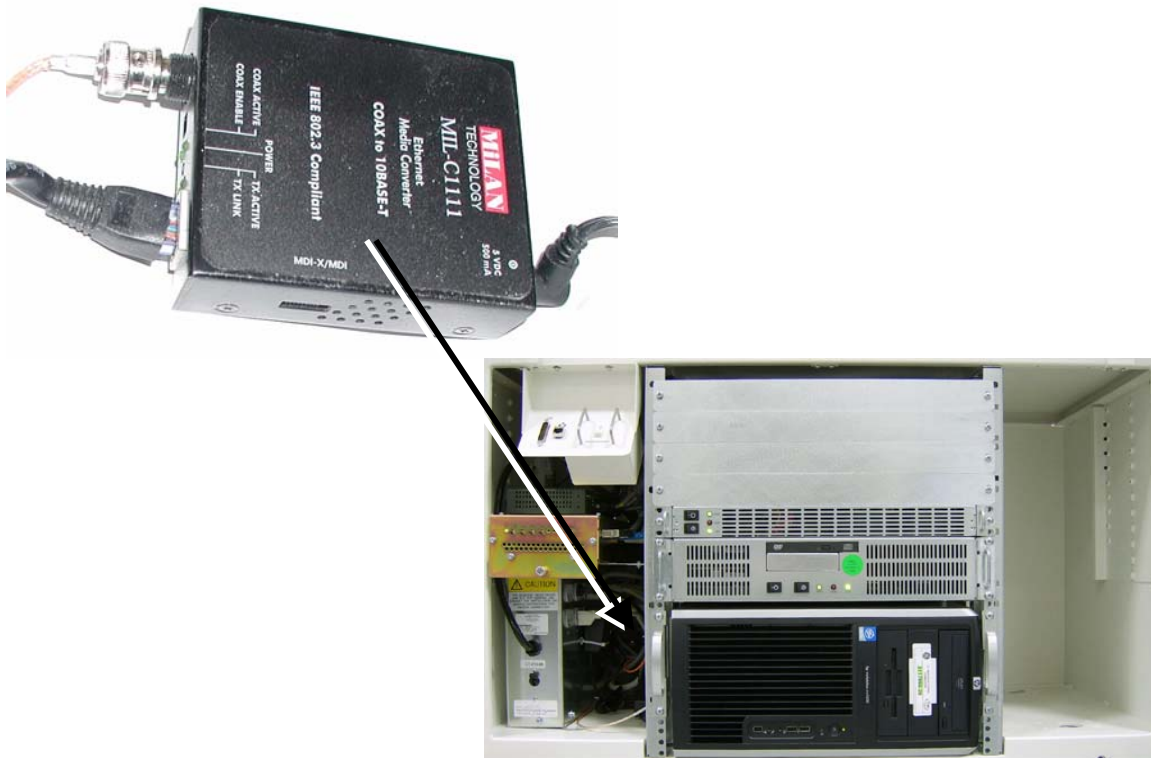


Illustration 2: Network Media Converter, Front Switch Setting

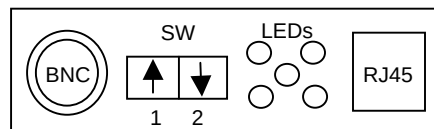

(Settings: SW 1—Up, SW 2—Down, no BNC terminator connected to BNC connector)

Illustration 3: Network Media Converter, Side Switch Setting

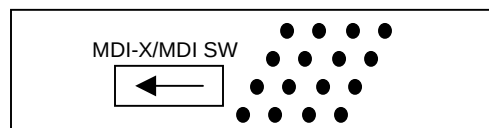

(Setting: “MDI-X” Position)

1.3 Daily Quality Reminder (Host Console Reboot)

Please advise your customer to reboot daily as part of their daily system quality procedures. It will prevent CT host console PC memory issues related to possible memory leakage.

2. Trackball PS/2 to USB Cable Converter

All PS/2 cable-based trackballs must be updated with the PS/2-to-USB cable converter included in this FMI kit, replacing the current black cable converter with this gray-colored version. Replacement is required for the PS/2-based trackball to work with the Linux OS 6.2.9 delivered by software collector 07MW11.10, distributed by this FMI.

If your trackball is PS/2-based, then replace its PS/2-to-USB cable converter as instructed below. If the trackball is USB-based, then skip this section.


NOTE: For PS/2 cable-based trackballs only, failure to replace the trackball's PS/2-to-USB cable converter may cause the CT host console to "hang" when loading, or after loading, 07MW11.10 software.

5. View the bulkhead connector at the top, rear of the host console. Refer to **Illustration 4** Host Computer Console Bulkhead for Trackball.

If the trackball connection is USB as in **Illustration 4** below, then skip this section and proceed to the next section; do not install the PS/2-to-USB Cable Converter.

If the trackball connection is PS/2, then replace the PS/2-to-USB converter cable perform next two steps.

Illustration 4: Host Computer Console Bulkhead for Trackball

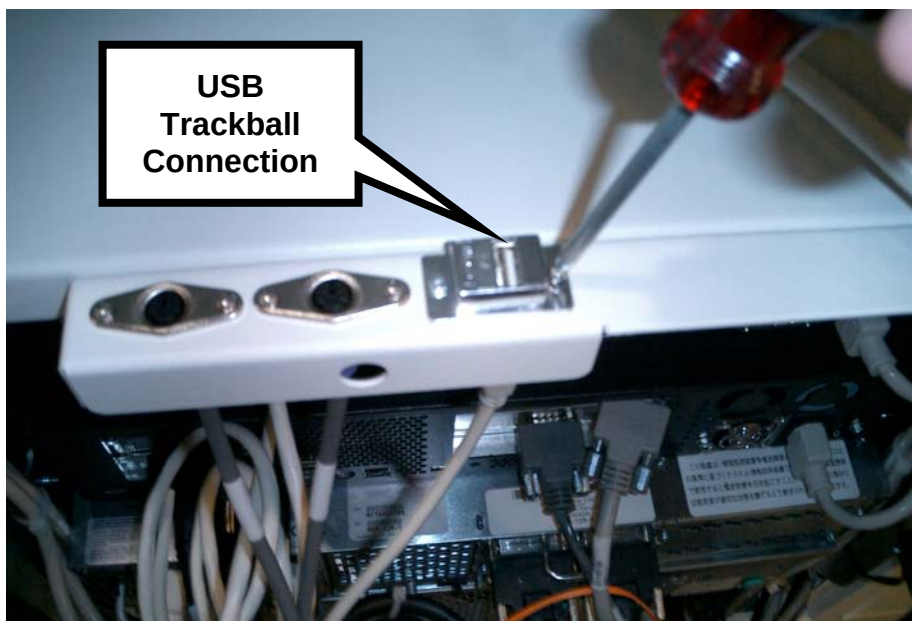

- a) Remove the rear cover to the host console to access the rear of the CT host computer console's HP PC.
- b) Replace the PS/2-to-USB cable converter (black) with p/n 5195935, the PS/2-to-USB cable converter (gray) that came in this FMI kit.

Refer to **Illustration 5** HP PC Rear View, PS/2-to-USB Cable Converter for an illustration of the PS/2-to-USB cable converter to be replaced.

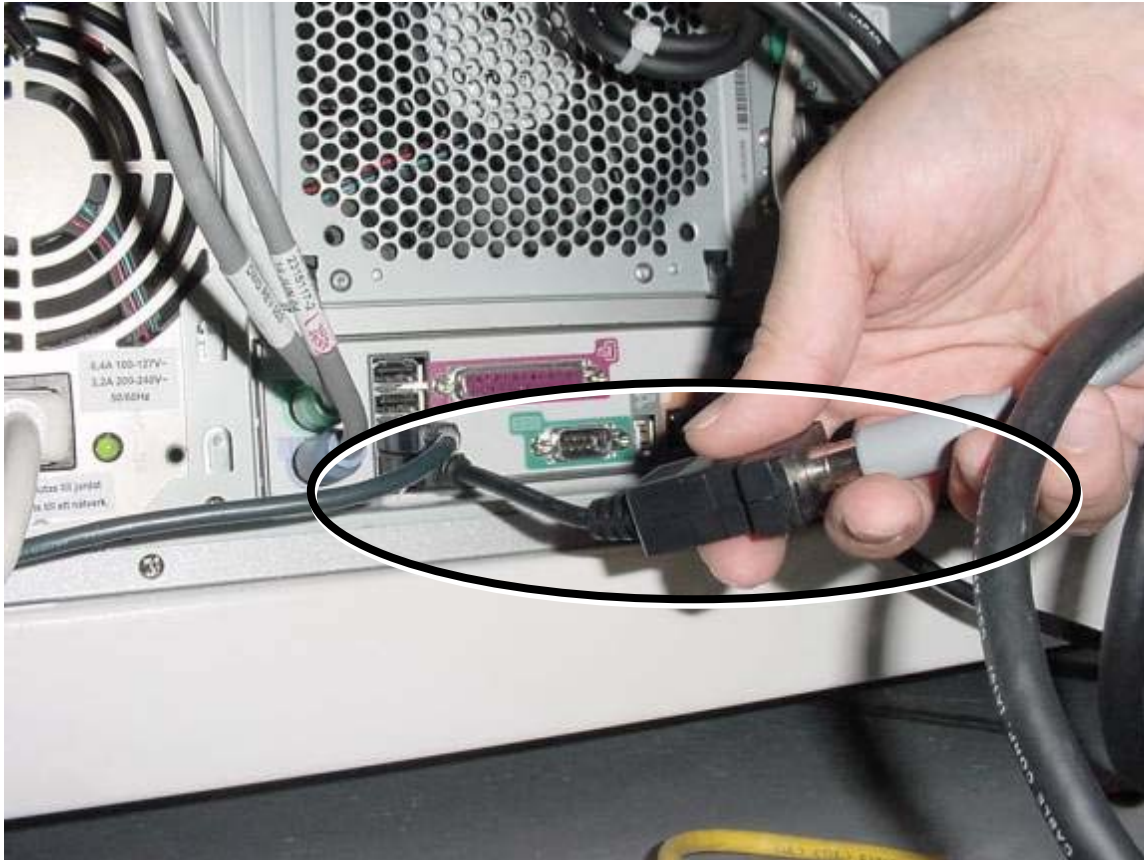


Illustration 5 HP PC Rear View, PS/2-to-USB Cable Converter

3. Software Load From Cold (LFC) Procedure

Software release 07MW11.10, collector p/n 5220443, will be loaded in this section.

07MW11.10 SOFTWARE ALREADY LOADED? (See "Prerequisites" section)

If the system is already running software release 07MW11.10, do not execute this FMI's load-from-cold procedure. Jump to Section 8, entitled "Service Pack 1.5 Software Load," to load Service Pack 1.5 software.

3.1. Need-To-Know Before Performing LFC

The following information is needed to correctly complete the 07MW11.10 load-from-cold (LFC) procedure:

1. In 07MW11.x LFC Procedure, Section 8.2 entitled "Check DARC Node Type," Step 6, only Westville DARC BIOS settings must be temporarily changed to support its software load. Jarrell DARCs do not require BIOS setting changes.
2. In 07MW11.x LFC Procedure, Section 8.2 entitled "Check DARC Node Type," DARC pictures are incorrect (IG is pictured). Refer to service information CDROM delivered by this FMI (p/n 5167277-200, Rev 6) → "Replacement" folder → "DARC2 Replacement" for correct picture of Westville and Jarrell DARC.
3. In 07MW11.x LFC Procedure, Section 8.4 entitled "OS / Applications Installation for DARC from Host Computer," Step 1.d, use the following command to always load the OS and applications onto the DARC:
`./load_darc_subsystem -restart`

If the "-restart" switch is not used in the command, the OS may not be loaded onto the DARC.

For more information, refer to service note SN07062211a, entitled "load_darc_subsystem -restart Script for DARC Software Reloads." It can be found at the following link:

http://www.ct.med.ge.com/tech_pubs/fct/system/servicenotes.htm

4. Per 07MW11.x LFC Procedure, Section 9 entitled "Service Software," the description listed may not be clear. It is at this point that Class C and Class M software should be loaded.
5. Per 07MW11.x LFC Procedure, Section 13 entitled "Install Software Options," Step 11, Table 6, software options must be reloaded in their proper sequence. eLicense-based software options reloaded/restored from the system's Save State DVD should be individually uninstalled, then individually reinstalled per the order listed in Table 6.
6. BIOS Floppies: BIOS load issues have occurred at a few early-installing sites. GEMS, nor the supplier, have been able to fully understand the root cause. Preventative action has been taken with the supplier, but recommend FEs be prepared by downloading BIOS 2.10 & BIOS 1.13 from Support Central to their laptops BEFORE beginning this LFC.

BIOS 2.10 and BIOS 1.13 files may be retrieved from the Support Central / CT Online Support website at

http://supportcentral.ge.com/products/sup_products.asp?prod_id=15768, then clicking *Tools / Patch Downloads / HP Computer Bios*

7. Autovoice loudness and feedback issues (squeal / noise) may occur during, or after, the load of 07MW11.10 software. If you experience these Autovoice issues, refer to Appendix C entitled "Autovoice Set Up: Volume and Feedback Issues" for corrective action.
8. Service Pack 1.5 software (p/n 5234957-15 inside software collector p/n 5234959-15 shipped in this FMI kit) must be installed for VoxTool / Reformat screens to be properly translated into the language for which the system has been configured. If not installed, the operator will see a mix of English and the selected language in VoxTool / Reformat screens.

When Service Pack 1.5 is installed, all previous software updates will be loaded (such as Service Pack 1.0). Refer to section 8 entitled "Service Pack 1.5 Software Load" of this FMI instruction.

3.2 Perform LFC

Perform the 07MW11.x software load from cold (LFC) procedure found on the "FX/Xtream Software and Upgrade Procedures" service information CDROM, p/n 5167277-200, Rev 6, shipped inside this FMI kit.

NOTE – Must Use 07MW11.x LFC: Failure to follow the 07MW11.x LFC procedure may cause increased customer downtime. The LFC procedure has been changed.

1. Insert the "FX/Xtream Software and Upgrade Procedures" service information CDROM into your FE laptop's DVD/CDROM drive.
2. Click on the "Software" folder (top-left side of display).
3. Click on "07MW11.x SW LFC" & complete the LFC procedure.

***** RETURN HERE WHEN THE LFC HAS BEEN COMPLETED. *****

4. Load Class C and Class M Software

To ensure proper operation of the InSite™ (includes InSite Update/Software Download) software, please load Class C and M in the following sequence:

- Load 07MW11.10 Class C software, CDROM p/n 5220446.
- Load 07MW11.10 Class M software, CDROM p/n 5220446-2.
- Re-initialize InSite™ Connectivity

4.1 Load Class C Software

Load Class C software per the instructions listed below.

If Class C and Class M software have already been loaded as part of the LFC procedure, then proceed to Section 5 entitled "Software Download / InSite Update File Ownership."

4. Open a shell.
5. Enter root as superuser:
 - a. At the ctuser prompt, type: **su** – [Enter]
 - b. At the Password prompt, type: **#bigguy** [Enter]
6. Insert the Class C software CD into HP computer DVD Drive. Wait for the disk read to STOP.
7. Install the Class C software on the system. At the root prompt, type: **install.classC** [Enter]
8. Close window after install.
9. Verify Class C is properly loaded on the system.
 - a. Open a Unix Shell
 - b. At the root prompt, type: **cd /var/adm** [Enter]
 - c. At the root prompt, type: **ls** [Enter]
 - d. Find the install.classc.log file. If file found, this means the Class C software load completed.
 - e. Remove the Class C software CD from HP computer DVD Drive.

4.2 Load Class M Software

If your system does not have InSite™ entitlement, do not execute this section and jump to Section 7 entitled “Save System State.”

If your system does have InSite™ entitlement, then load 07MW11.10 Class M software as follows:

10. Open a shell.
11. Enter root as superuser:
 - a. At the ctuser prompt, type: **su** – [Enter]
 - b. At the Password prompt, type: **#bigguy** [Enter]
12. Insert the Class M software CD into the HP computer DVD drive. Wait for the disk read to STOP.
13. Install the Class M software on the system. At the root prompt, type:
install.classM [Enter]
14. Close window after install
15. Verify Class M is properly loaded on the system.
 - a. Open a Unix Shell
 - b. At the root prompt, type: **cd /var/adm** [Enter]
 - c. At the root prompt, type: **ls** [Enter]
 - d. Find the install.classM.log file. If file found, this means the Class M software load completed.
 - e. Remove the Class M software CD from HP computer DVD Drive.

5 Software Download / InSite Update File Ownership

To enable usage of Software Download / InSite Update buttons in Service Desktop → Utility → Update Management, the ownership of three files listed below must be set to “ctuser.”

16. Open unix shell and verify that you are “ctuser” (not as “root”).
17. Type **touch /tmp/cl** [Enter]
18. Type **touch /tmp/cl1** [Enter]
19. Type **touch /tmp/cl2** [Enter]
20. Verify file ownership for cl, cl1, and cl2 is now “ctuser,”
 - a. Type **cd /tmp** [Enter]

- b. Type **ls -al cl*** [Enter]

Output should look similar to below, note the ownership is “ctuser.”

```
-rw-rw-rw- 1 ctuser users 9 May 31 11:38 cl  
-rw-rw-rw- 1 ctuser users 2 May 31 11:38 cl1  
-rw-rw-rw- 1 ctuser users 27 May 31 11:38 cl2
```

21. Close unix shell.

6 Re-Initializing InSite™ Connectivity

6.1 Restore System State To Retrieve Class M Configuration Parameters

Restore 'Reconfig' part of the Restore System State to recover all previous InSite™ and Class M configuration parameters.

6.2 Enable InSite™ Interactive Platform Configuration Utility

Enable the InSite™ Interactive Platform Configuration utility as follows:

22. Open a shell.

23. Enter root as superuser:

- a. At the ctuser prompt, type: **su -** [Enter]
- b. At the Password prompt, type: **#bigguy** [Enter]

24. Launch the IIP Admin Configuration utility. At the root prompt, type:
iipadmin config [Enter]

The system displays the InSite™ Interactive Platform Configuration utility with the Proprietary Agreement Acceptance Screen.

25. Click on **ACCEPT** three times to display the Information Screen.

26. Click on the **ProDiags** Tab to display the ProDiags Screen.

Click on the **Custom** button to launch the independent ProDiags

This will launch ProDiag Utility.

27. There is no change to ProDiags. However, it is recommended that you review ProDiags Schedule and update, if needed, according to your particular site's normal operation. The Schedule Tab is used to modify Scheduled Tasks.

As an example, use the following sequence to update the scheduled time (hour and minute) of the **idm** and **idm_send_data** tasks

1. Click on **Hour** and **Minute** fields of **idm_send_data** task and select an hour and minute from the pull down lists between 9 AM and 5 PM.

After each selection, click **Yes** on the "Data has been modified. Please confirm it" pop up.

2. Click on **Hour** and **Minute** fields of **idm** task and select an hour and minute from the pull down lists to schedule 5 minutes before the **idm_send_data** task.

After each selection, click **Yes** on the "Data has been modified. Please confirm it" pop up.

3. When **idm** and **idm_send_data** Hour and Minute selections are complete, click on **Apply** to finalize the schedule updates.

Click **Yes** on the "Are you sure you want to apply changes?" pop up.

Task	Iterations	Day	Hour	Minute	Type	Time Slot
PassChg	1	Daily	05 AM	30	Background	
pd_HouseK...	1	Daily	04 AM	30	Background	
pd_ASC_No...	1	Daily	Hourly	28	Background	
healthng	1	Daily	04 AM	45	Background	
IIP_HouseK...	1	Daily	06 AM	12	Background	
IIP_HouseK...	1	Daily	12 PM	12	Background	
IIP_HouseK...	1	Daily	06 PM	12	Background	
IIP_HouseK...	1	Daily	12 AM	12	Background	
idm_send...	1	Daily	07 AM	00	Background	
idm	1	Daily	Hourly	00	Background	

Buttons: New, Delete, Reload, Apply, Help

Status: Reading data... Done

28. Close the ProDiags utility browser by selecting **Close** from the upper left browser window pull down menu. The ProDiags browser should close and the InSite™ Interactive Platform Configuration utility should display the HealthPage Tab as shown in the figure on page 19.

6.3 Configure HealthPage


After completing the ProDiags setup, the IIP Admin Configuration utility automatically switches to the HealthPage.

The Restore State process (after loading Class M software) restored the following information for the HealthPage Reports:

- All e-mail addresses of targeted recipients that were previously loaded
- HealthPage report contents that was previously configured
- HealthPage report schedule that was previously configured

Unless you have recent changes to any of these three items, you can skip the HealthPage Tab of the IIP Admin Configuration utility and select the Device Connection Tab. If you have changes that you want to make (e-mail addresses,

report content, report schedule), follow the instructions on the IIP Admin Configuration utility.

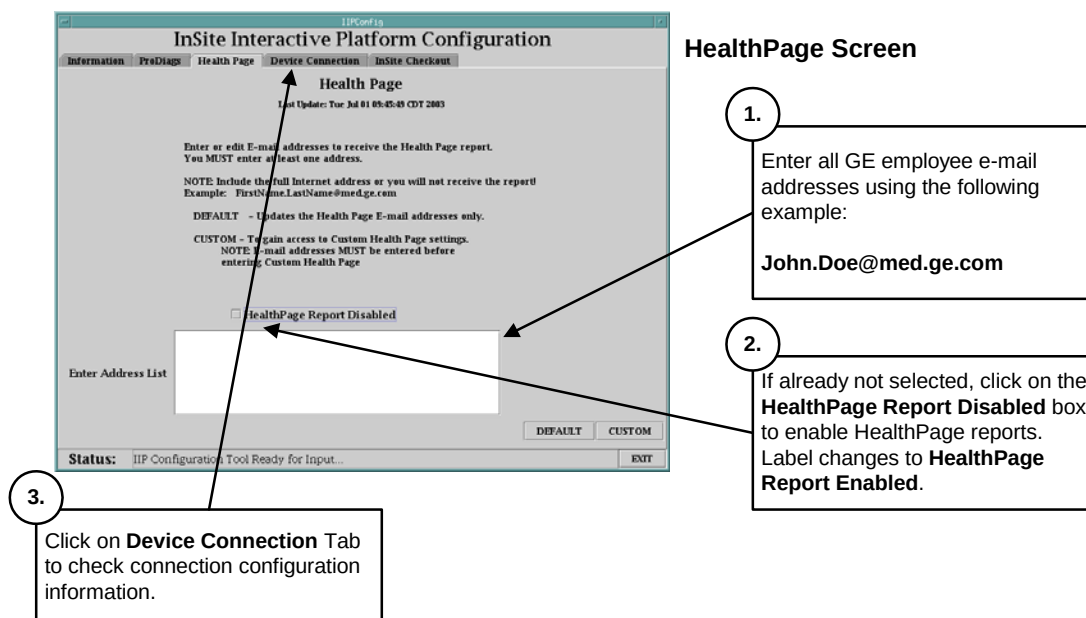

6.4 Check Device Connection Setup

When selecting the Device Connection Tab, the IIP Admin Configuration utility automatically displays the Modem Setup Screen.

Make sure that the modem set-up parameters are correct. If not, change the parameters to the correct ones. See Appendix A.

6.5 Check BroadBand Connection Setup

If your system was configured for BroadBand connectivity, the Restore State process (after loading Class M software) restored all the BroadBand setup parameters. If not, enter the proper IP address of the InSite™ Gateway. See Appendix B for instructions.

6.6 Completing the InSite™ Re-Initialization

Follow the instructions provided below to exit the IIP Admin Configuration Utility.

Insite Checkout Screen

1. Click on the Exit button of the IIP Admin Config utility to exit the utility.

2. Click on OK on the Exit Confirmation pop-up to exit the IIP Admin Config utility and run the ConfigLink utility.

3. After the shell displays the root prompt, reboot the system to load Insite connectivity parameters. As the system is coming up, you should see the iLinq icon (if the original Checkout included iLinq) as a selectable desktop on the right console monitor.

InSite Interactive Platform Configuration - InSite Checkout

Last Update: Tue Jul 01 09:45:59 CDT 2003

To Perform an InSite Checkout:
1. Verify all Tab information has been updated.
2. Select one of the buttons below.

CHECKOUT NOW - Use this to set system ready for OLC to complete the checkout process.
Requires the FE to call support center to do checkout.

DISABLE INSITE - InSite Disabled until IIP Configuration is run again and the Device Connection is re-configured.

AUTO CHECKOUT - Use this if InSite for this site has been checked-out successfully at least once in the past and PreDiags has been configured.

Host List			
127.0.0.1	localhost.localdomain	localhost	
1231.15.62	iblah2		
1929.220.1	c099 ac CT99 setb	CT99_OCB ac	
1929.220.11	CT99_ice ice ICE		
1929.220.12	ice-sim-gate		
1929.220.13	ice-pant		
1929.220.14	ice-signalb		
1929.220.1	shreerig		

Buttons: CHECKOUT NOW, DISABLE INSITE, AUTO CHECKOUT, EXIT

Status: Reading Modem Configuration Files...

Exit Confirmation

Configuration Last Updated Summary

Proactive Diagnostics	: Mon Oct 06 18:27:02 CDT 2003
Health Page Email	: Thu Oct 09 14:55:20 CDT 2003
Health Page Options	: Tue Oct 07 14:33:36 CDT 2003
Health Page Dates	: Tue Oct 07 14:33:36 CDT 2003
Health Page Reporting	: ENABLED
Connection Configuration	: Thu Oct 09 15:10:23 CDT 2003
Checkout Status	: Wed Oct 08 10:15:57 CDT 2003

Press OK to Exit OR Press Cancel to abort exit

Buttons: OK, Cancel

6.7 Run ConfigLink

Run the InSite™ **ConfigLink** script to populate all configuration information required to fully restore connectivity parameters on the system.

29. Open a shell.

30. Enter root as superuser:

- a. At the ctuser prompt, type: **su -** [Enter]
- b. At the Password prompt, type: **#bigguy** [Enter]

31. Change directory to InSite™. At the root prompt, type:

cd /usr/g/insite [Enter]

32. Run the ConfigLink script. At the root prompt, type:

ConfigLink [Enter]

Wait for the shell to return the root prompt to ensure the **ConfigLink** process completes. When complete, exit root and close the shell.

6.8 Verify telnet in PNF (Product Network Filters) is Off


The default setting in PNF allows telnet from all IP addresses. To make the system more secure, turn off telnet via the commands below. If telnet is required for systems within the hospital network, individual IP addresses can be added as needed.

33. Open Common Service Desktop

34. Click on Configuration Tab

35. Click on Configure PNF.

36. Follow steps 1-4 shown in picture on the next page to turn off Product Network Filters (PNF).

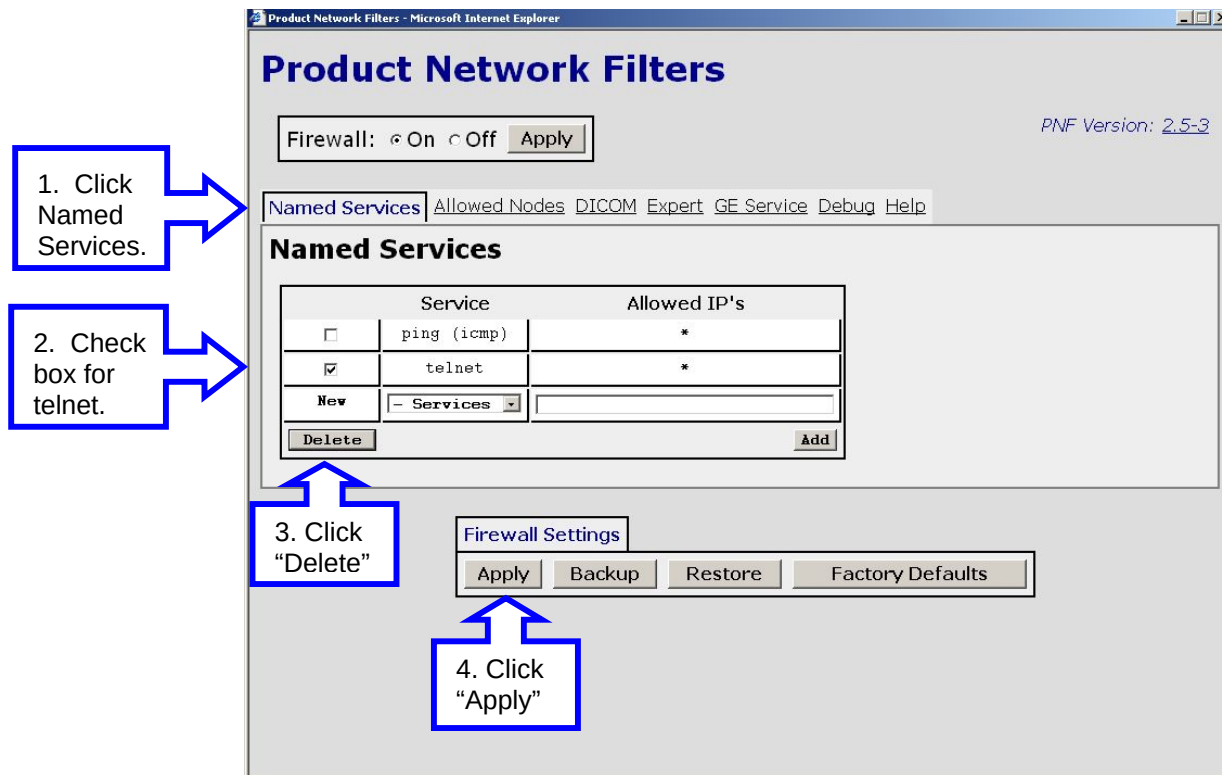

6.9 Reboot System

37. Click on the pink Shutdown button.

38. From the Shutdown selection pop-up, click on RESTART, then OK, to reboot the system.

7. Save System State

39. Save Reconfig portion of the system state using the Save System State utility.

8. Service Pack 1.5 Software Load

Service Pack 1.5 CDROM (p/n 5234957-15 inside software collector p/n 5234959-15 shipped in this FMI kit) must be installed for VoxTool / Reformat screens to be properly translated into the language for which the system has been configured (only applies to non-English language-configured systems).

Check if the Service Pack 1.5 software has been loaded when the 07MW11.10 LFC procedure was executed (Section 24 entitled "Finalization," Step 15 of the LFC procedure). Load it now if it was not loaded during the LFC.

40. Check if Service Pack 1.5 software is loaded.

a) Open Unix shell

- b) Login: **ctuser** [ENTER]
- c) Type: **showprods | grep -i servicepack** [ENTER]

Example of returned string:

07MW11.10_07MW18.4_07BW08.4_GEHCSERVICEPACK1.5
08/20/2007 07MW11.10/07MW18.4/07BW08.4 Service Pack 1.5 ver 1
rel 1.2

41. If the returned string does include "Service Pack 1.5," then jump to Section 9 entitled "System Checkout" because the service pack software is already installed.

If the returned string does not include "Service Pack 1.5," then install Service Pack 1.5 software now. Install it as instructed in Step 41a below.

NOTE: Service pack installation will fail if an earlier service pack is currently installed.
Uninstall the current service pack now, if one is viewed when performing the above **showprods** command (example, "07MW11.10_servicepack 1 ver 1 rel 1.1"). Refer to the service pack CDROM label for uninstall instructions.

- a) Retrieve Service Pack 1.5 CDROM, p/n 5234957-15 inside software collector p/n 5234959-15, from the FMI kit
- b) Insert the CDROM into the HP PC DVD drive
- c) Shutdown Apps
- d) Open Unix shell
- e) Login: **su -** [ENTER]
- f) Type: **patch_install -c** [ENTER]
- g) Type: **reboot** [ENTER]

9. System Checkout

42. Perform the following sanity scan to check system functionality (Abdominal or Head) -- 50 images at 0.5 secs rotation at 120kV, 40mA

The system should scan and image without any errors

10. Rear Console Cover

43. Re-attach the rear console cover.

11.Finishing Up

44. Retrieve the correct set of documents, from this FMI kit, to be left on site with the imaging department leader.

Use the tables below to determine which documents and CDROMs are specific to your system. After you have located the documents specific to your system, leave them on site and discard all other remaining documents.

LightSpeed 1X	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide – English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

LightSpeed 2X	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide – English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

LightSpeed 3X (Includes HiSpeed QX/i)	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide – English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

LightSpeed 4X	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide -- English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

LightSpeed 5X 100kW / 80 kW	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide -- English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

LightSpeed 5X HP60 4/8 slice	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide -- English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

LightSpeed 5X HP60 16 slice	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide -- English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

LightSpeed 5X RT 4 slice	
2369740-247	CT LIGHTSPEED SERIES CD
5119787-100	LS Series Quick Step Guide -- English
5119787-247	LS QuickStep Guide
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

BrightSpeed 4 slice Edge / 8 slice Excel	
5141172-247	LRG Operator Manual
5159870-100	BrightSpeed QuickStep_English
5159870-247	BrightSpeed QuickStep Guide - Multi-language CD
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

BrightSpeed 16 slice Elite	
5141172-247	LRG Operator Manual
5159870-100	BrightSpeed QuickStep_English
5159870-247	BrightSpeed QuickStep Guide - Multi-language CD
5198350-100	2007 Common Platform Software Release Application Tips & Workarounds, Paper Version
5198350-247	2007 Common Platform Software Release Application Tips & Workarounds, MultiLanguage CD
5231794-148	Product Matrix Manual

45. If the system is installed in an English-speaking country, retrieve and leave at site with the imaging department leader the following additional application documents:
- a. 5183783-247for all system types
 - b. 5252415-199for all BrightSpeed scanners only

46. If the system is installed in a Japanese-speaking country, retrieve and leave at site with the imaging department leader the following additional application documents:
- a. 5183783-247for all system types
 - b. 5198350-140for all system types
 - c. 5119787-140for all LightSpeed scanners only
 - d. 5159870-140for all BrightSpeed scanners only
 - e. 5252415-199for all BrightSpeed scanners only
47. If system is installed in a Chinese-speaking country, retrieve and leave at site with the imaging department leader the following additional application documents:
- a. 5183783-247for all system types
 - b. 5198350-141for all system types
 - c. 5119787-141for all LightSpeed scanners only
 - d. 5159870-141for all BrightSpeed scanners only
 - e. 5166779-141for all BrightSpeed scanners only
 - f. 5252415-199for all BrightSpeed scanners only


12.FMI Closure

48. Deliver the appropriate application documents to the imaging department leader to support the use of the scanner running software release 07MW11.10.
49. Remove FMI Instructions from the customer site.
50. Close FMI per regional procedure and close.

APPENDIX A

Modem Setup

When selecting the Device Connection Tab, the IIP Admin Configuration utility automatically displays the Modem Setup Screen as shown below.

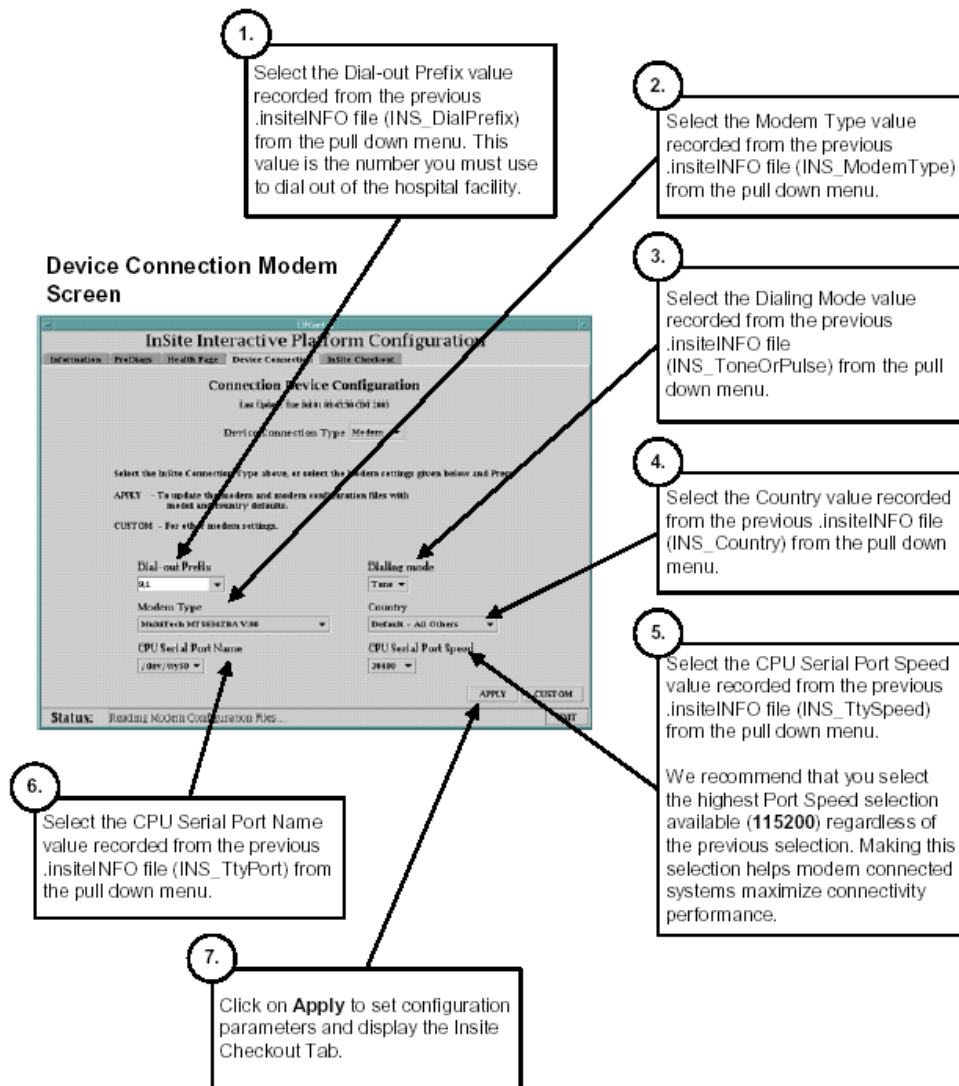

If your system was configured for Modem connectivity, the Restore State process (after loading Class M software) restored all the Modem setup parameters. Verify with the InSite™ information that you recorded before starting the loading process.

- If there were no Modem parameters in the .InSite™INFO file but there were BroadBand parameters, skip to the Check BroadBand Device Connection Setup procedure.
- If there were Modem parameters in the .InSite™INFO file but the parameters displayed in the Modem Setup Screen do not match what you previously recorded, enter the correct values as required.

NOTE: The Modem configuration parameters should be pre-loaded from the Restore State process. Do not change settings unless these values were not properly restored.

Appendix B

Broadband Connection Set Up

If your system was previously BroadBand InSite™ connected, click on the Network button of the IIP Admin Configuration utility Device Connection Tab to display the Network Setup Screen.

If your system was configured for BroadBand connectivity, the Restore State process (after loading Class M software) restored all the BroadBand setup parameters. If not, enter the proper IP address of the InSite™ Gateway as shown in the figure below.

NOTE: The BroadBand Network configuration parameters should be pre-loaded from the Restore State process. However, the on-site FE who performed the original InSite™ Checkout of the system may have used the hospital's imaging LAN default gateway IP Address (selection was available in the previous release of IIP).

DO NOT USE THE HOSPITAL WORKING LAN DEFAULT GATEWAY IP ADDRESS AS THE InSite™ GATEWAY IP ADDRESS OTHER THAN AS A TEMPORARY MEANS TO CONNECT TO THE ON LINE CENTER!! Contact the hospital network administrator to determine the right gateway address.

Device Connection Network Screen

1. Enter the IP Address of the Insite Gateway to GE Online Center as recorded from the .insiteINFO file (INS_GatewayIP). Do not enter leading "zero's" into the four fields.
2. Click on **Apply** to have the IIP Admin Config utility ping the entered IP Address, load the network parameter, and display the Insite Checkout Tab.
3. Click on **OK** on Checkout pop-up


Appendix C

Autovoice Set Up: Volume and Feedback Issues

Autovoice low-volume and feedback (noise) issues may be resolved by changing two audio mixer settings via the “alsamixer” graphical user interface (GUI). If the system exhibits any of these issues, follow the instructions listed below in the given order:

1. If Autovoice volume is too low, select “Autovoice Volume” tool (“audioSetting” tool) from the Tool Chest menu to increase it as needed, and click the “Save,” then “OK” buttons to save to file.
2. If Autovoice feedback issues exist (increasing the volume setting may create feedback), then
 - a. Open a unix shell.
 - b. Type: **alsamixer** [Enter]

The following screen will appear:

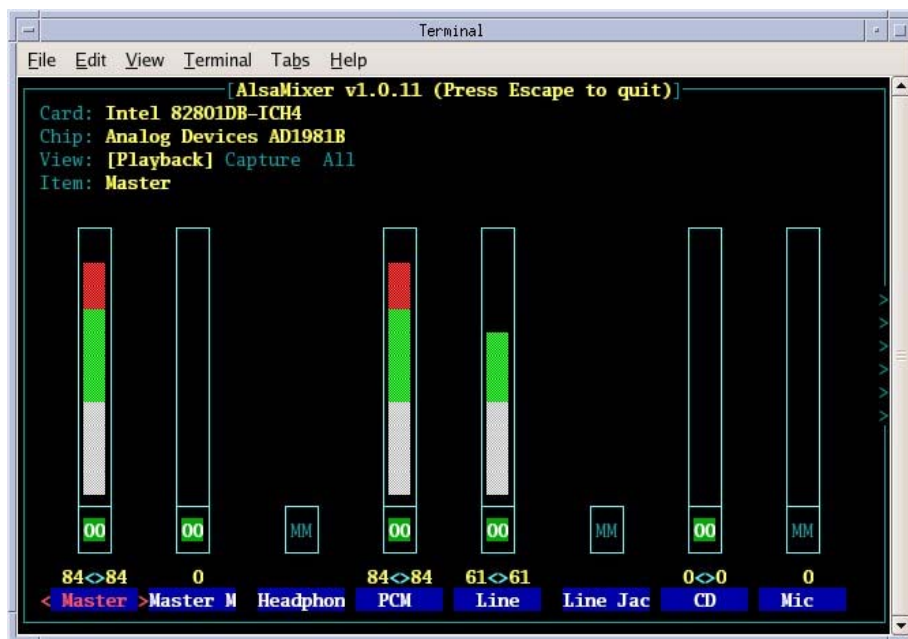

NOTE: Settings shown here may not depict the final settings for your console. Settings vary from console to console.

- c. Press the right-arrow key → to move cursor selection to the **PCM** slider.
- d. Press the up-arrow key ↑ to increase the **PCM** setting to 100%.
- e. If feedback is heard, press the right-arrow key → to move cursor selection to the **Line** slider.

Press the down arrow key ↓ to decrease **Line** setting until feedback is no longer heard.

Do not mute **Line** device; it is needed for Autovoice Record.

f. When completed, press the Esc key to exit AlsaMixer.

g. If the **Line** setting was changed, then STOP → you have completed Autovoice setup.

If the **Line** setting was not changed, then select “Autovoice Volume” tool from Tool Chest menu and click the “Save,” then “OK” buttons to save the **PCM** settings to file.

NOTE: If the **Line** setting was changed and then Autovoice Volume settings saved, the **Line** value set via AlsaMixer was overwritten by the “Save” and set back to “61.” Steps 2e and 2f must be repeated to adjust the **Line** setting.