

**ADMINISTRATIVE INFORMATION****FMI No: 25402****FMI Title: VCT VIG EMC Update**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

- ☒ **Staged FMI** **Customer List** provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
- ☒ **Parts** and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
- ☐ **Parts Ordering:** GEMS - Americas, Contact your Region logistics coordinator
GEMS - Europe, Order FMI kits from GPO/EPL Buc
GEMS - Asia, Contact your region logistics coordinator
- ☐ **No Parts** are required to complete this FMI.
- ☐ **Batch with** _____
- To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
- ☐ **Non-Staged** **No Customer List is provided.** Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS – Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:

This FMI applies LightSpeed 7.x VCT and Discovery VCT using Jarrell VIG model# 5159834-2. If the VIG model number is 5159834-3, this FMI does not need to be applied. This procedure must be done on all VIG's in the console on site.

Unit Tracked: LightSpeed 7.x and Discovery VCT systems**FMI Instructions Part No.: 5231915-1EN****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS – Europe	
GE Medical Systems Insite Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

FMI No.: 25402

FMI Title: VCT VIG EMC Update

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments

Table of Contents

Table of Contents	1
Revision History.....	2
Effectivity	2
Purpose	2
Furnished Materials.....	4
Tool Required	4
Personnel Requirements.....	4
Procedure	5
POWER-OFF (SHUT DOWN) THE CONSOLE.....	5
REMOVE VIG NODE FROM CONSOLE	5
VIG VRAC REWORK	6
INSTALL THE VIG NODE INTO CONSOLE	25
POWER-ON THE OPERATOR CONSOLE	26
VERIFY V/DARC NODE IS UP.....	27
VERIFY REPLACEMENT VIG NODE	27
VERIFY RECON DATA PATH DIAGNOSTIC	29
Finishing Up.....	29
Appendix A – ethtool T/S Note	30
Appendix B – Verifying Jarrell VIG Hardware	31
HARDWARE VERIFICATION	31
REMOTE VERIFICATION.....	32

Revision History

Revision	Date	Reason for Change
1	6/26/2007	Initial release.
2	9/5/2007	Updated wording in effectivity. Added ESD Note in VIG VRAC Rework Procedure.
3	10/18/2007	Added Appendix B for determining Jarrell vs. Westville VIG

Effectivity

This FMI affects all LightSpeed 7.x VCT (64, 32, Select) shipped with Jarrell VIGs (5159834-2), prior to 04/27/2007. This FMI also affects all LightSpeed 7.x systems that have ordered the Jarrell VIG (5159834-2) through FEMC. All Jarrell VIG's with model# 5159834-2 require rework. The model numbers tracked in this FMI are:

System	Model#	Description
LightSpeed VCT64	5124069	VCT Gantry
	5124069-2	VCT Gantry
	5124069-3	VCT Gantry
	5124069-4	VCT Gantry
	5124069-5	VCT Gantry
	5124069-6	VCT Gantry
LightSpeed VCT32	5126093	VCT PRO32 Gantry
Discovery VCT	5158924	DSTE VCT Gantry
	5178168	DRX VCT Gantry
	5193989	DRX VCT Gantry

For detailed instructions to check if you have a Jarrell VIG, see *Appendix B – Verifying Jarrell VIG*.

Purpose

During an internal GEHC manufacturing audit, it was found that the electromagnetic radiation emissions, from the console, were +2.5 dB over the legal limits defined in CISPR 11. Investigation determined that the rear of the console, more specifically the VRAC2 card installed in the IG enclosure, had radiated emissions above the IEC 60601-1-2 requirements.

Further investigation has determined the cause of excessive radiation emission is linked to improper grounding of the VRAC circuit board and the chassis design that contains the VRAC2. A new mounting scheme and faceplate has been developed by the engineering team to address

both issues. The redesigned and modified Jarrell VIG assembly will be assigned the new part number 5159834-3.

Furnished Materials

Qty	Description	Part #
1	FMI Document	5231915-1EN
3	Rating Plate Label	46-302200P12
1	Sharpie Fine-Point Perm. Marker	2168609
1	VIG Rework Kit	5248196
1	3/8" Hex Nut Driver	MSP-3/8NUTDRIVER
1	Anti-Static Bag 8X12	MSS-ESDBAG-8X12
1	Anti-Static Bag 6X8	MSS-ESDBAG-6X8
13	Silver Lid Screw	MSC-M3/.5/4/Z
3	VRA2C Mounting Bracket	FAB-BRKT0569
3	PCI Slot Cover Bracket	FAB-COVERPCI0570
1	VRAC2 Installation Jig	FXT-JIGVRAC
7	VRAC2 Mounting Bracket Screw	MSC-6/32/.1875/Z
10	PCI Slot Cover Bracket Screw	MSC-M3/.5/5MM/B
10	VRAC2 Screw	MSC-M3/.5/5MM/Z
7	3/8" Hex Standoff	MSC-STNDFF 3/8-.5
10	Spacer	MSP-SPACER03125

Tool Required

Qty	Description
1	Standard FE Toolkit
1	ESD Workstation Monitor
1	ESD Wrist Strap
1	ESD Anti-Static Mat
1	Flat Head Screw Driver

NOTE: ESD Workstation may not be applicable to all poles. Field personal are to follow the standard ESD practices and procedures for this FMI.

Personnel Requirements

"2" hours labor and "1" hours of travel

Procedure

NOTE: This procedure only shows the instructions to perform this FMI on one VIG. All VCT systems have the option of having up to three VIG's. This FMI must be performed on all VIG's in the console on site.

This FMI is only applies to Jarrell VIGs. For assistance identifying a Jarrell VIGs, see *Appendix B – Verifying Jarrell VIG*.



Warning: Static Sensitive Parts

This FMI requires field personal to handle circuit boards. The following Electrostatic Discharge must be used while performing this FMI:

- ESD Workstation Monitor and/or Wrist Strap
- Anti-Static mat

Power-Off (Shut Down) the Console

1. Select one of the following methods to Power OFF the Operator Console:
 - a. If Applications are up, click on the [Shut Down] icon and select [Shutdown].
 - b. If Applications are down, open a Unix Shell at the Toolchest. Type:
`{ctuser@hostname}: halt <ENTER>`

The Operator Console monitor will display a 'System halted' message when it is acceptable to power OFF the Operator Console.

2. Power **OFF** the Operator Console at the front panel switch.

Remove VIG Node from Console

NOTE: This FMI is only applies to Jarrell VIGs. For assistance identifying a Jarrell VIGs, see *Appendix B – Verifying Jarrell VIG*.

1. Remove the front and rear Operator Console covers.
2. Remove the rear Ethernet cable connection(s) to the VIG Node(s) being reworked. If more than one cable is being disconnected, verify each cable label is present and clearly marked; if necessary, add a label for clarity.
3. Remove the power cord at the rear panel of VIG Nodes.

4. Remove and retain the two (2) front mounting screws and slide out the VIG Node(s).

VIG VRAC Rework

NOTE: Work on the VIG must be performed with the VIG placed on the ESD mate, with proper grounding. Boards should be placed in ESD bags when removed from the VIG chassis!

1. Cut the Warranty Seal along the seam of the lid.




2. Remove and discard the 4 screws securing the lid, then remove the lid.




3. Remove and discard the 2 PCI bracket mounting screws and 4 washers using a P1 driver bit.


NOTE: These screws are secured with Loctite, but they will come off.

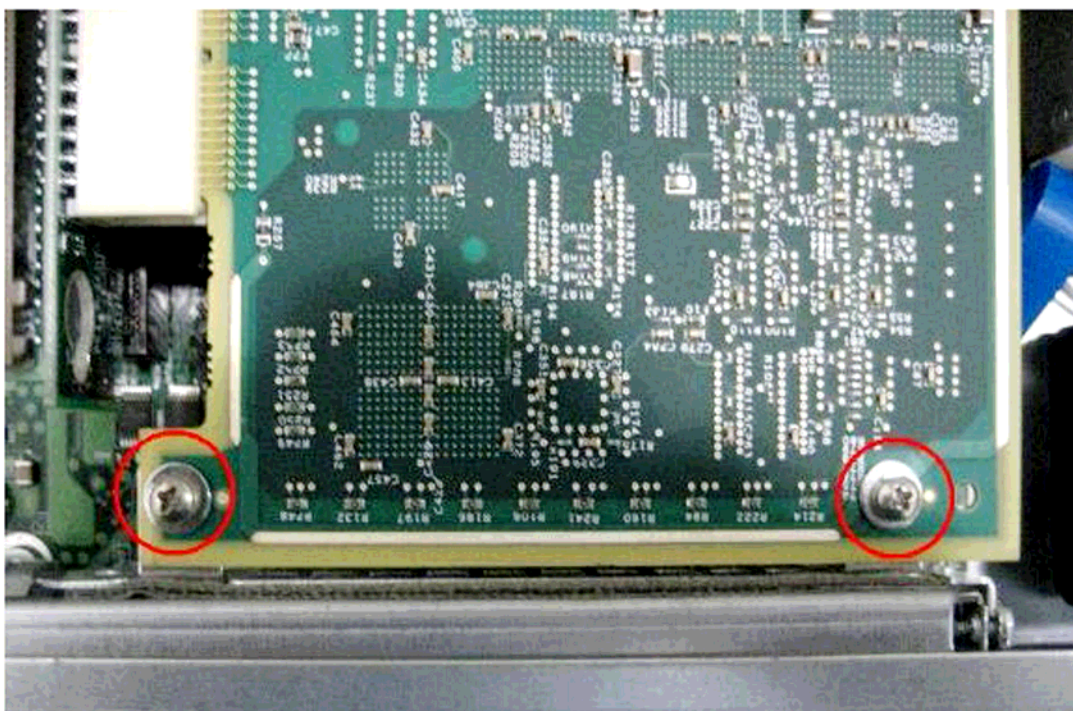


4. Remove and discard the screw from the corner of the VRAC card.

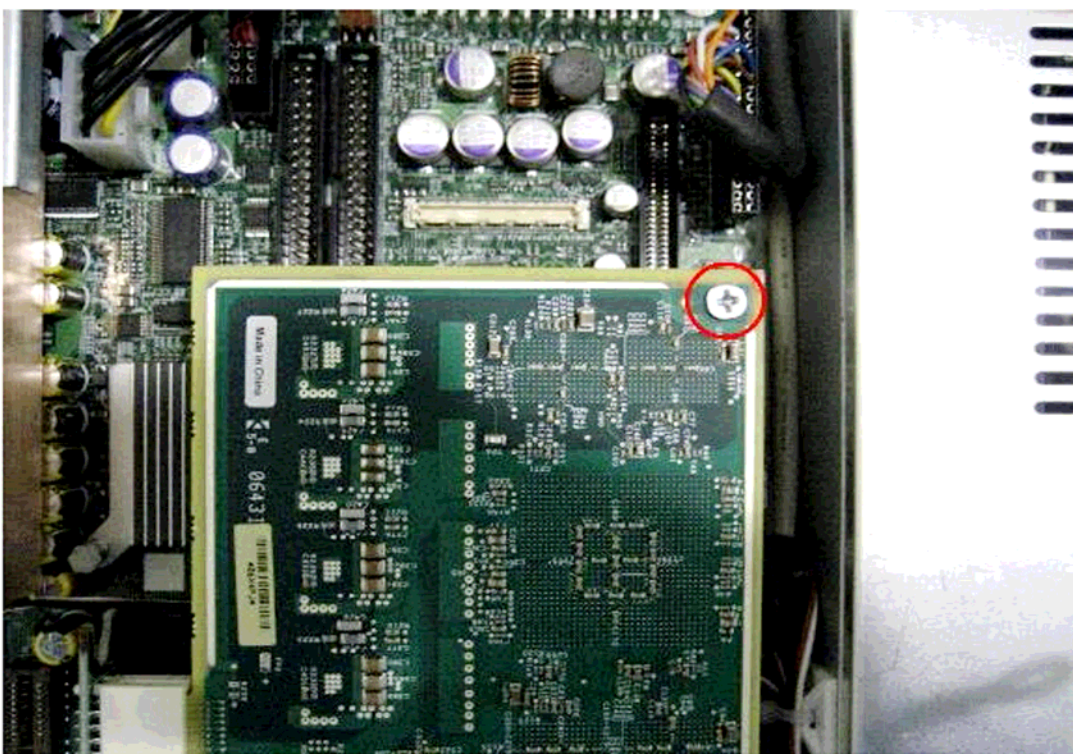



5. Remove the VRAC card and riser from the system by pulling straight up on both ends of the riser.

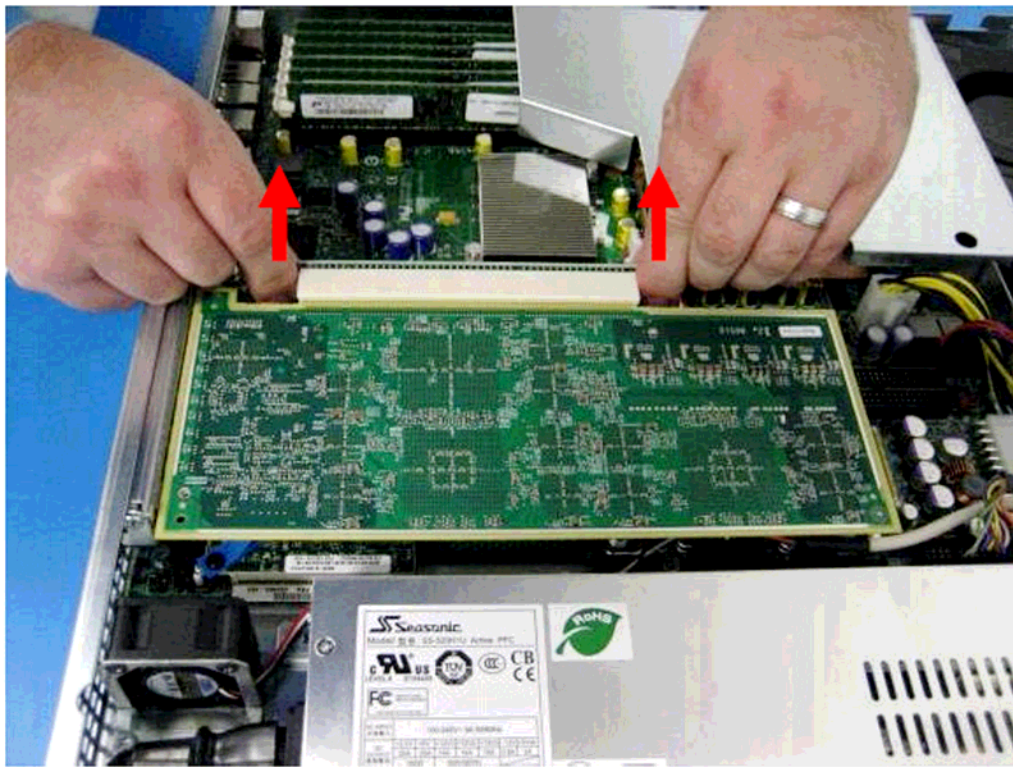


6. Carefully disconnect the VRAC card from the riser without touching the gold connectors.

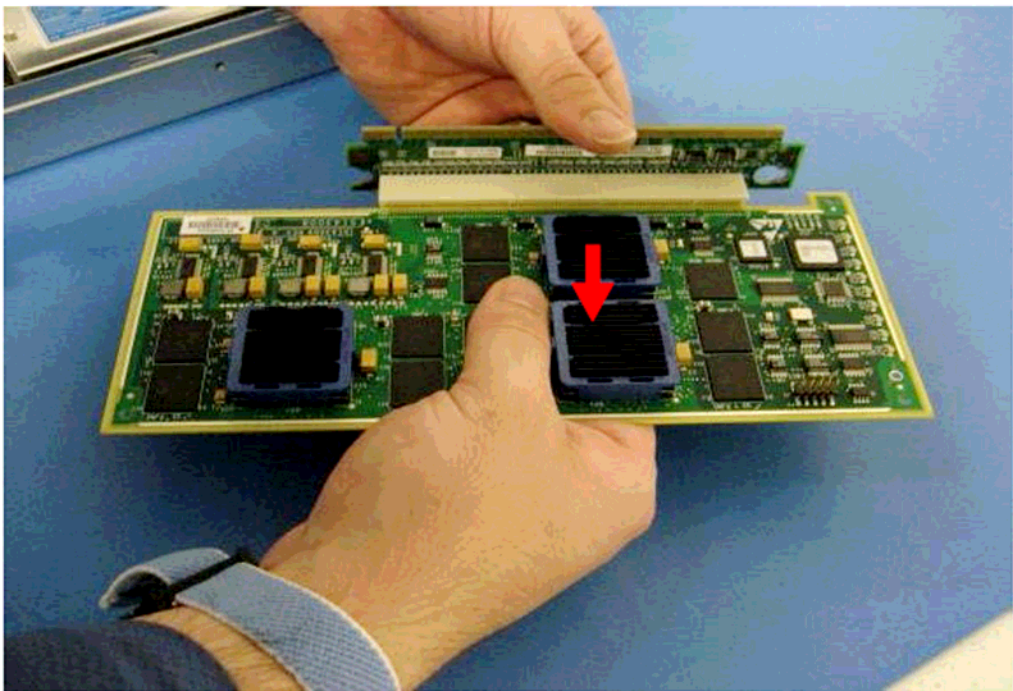

7. Place both components in the ESD-safe bags provided with the FMI (P/N .MSS-ESDBAG-8X12 and MSS-ESDBAG-6X8).

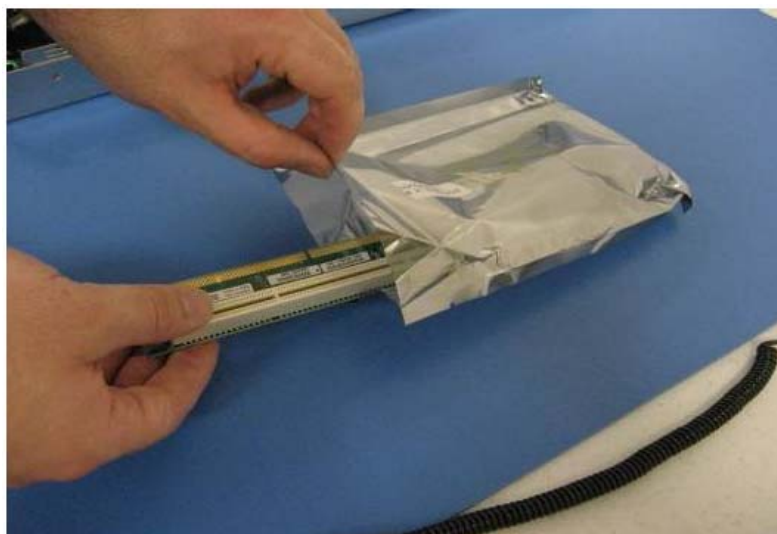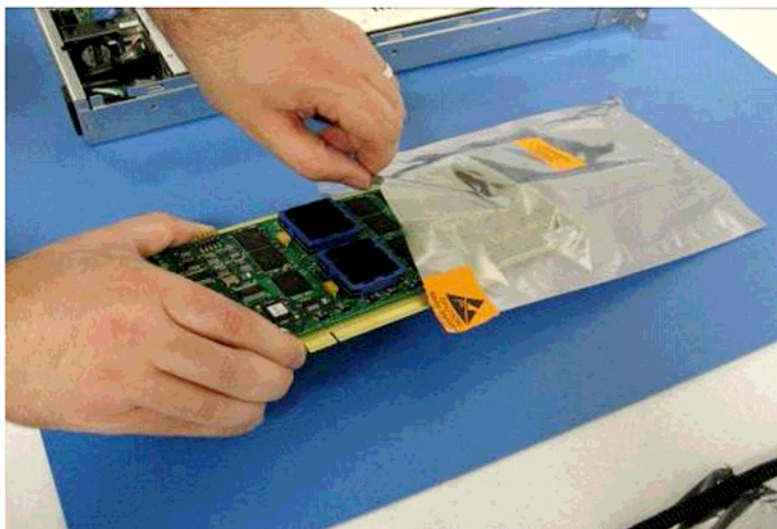

8. Remove and discard the 3 screws and PCI mounting bracket from the rear of the chassis.

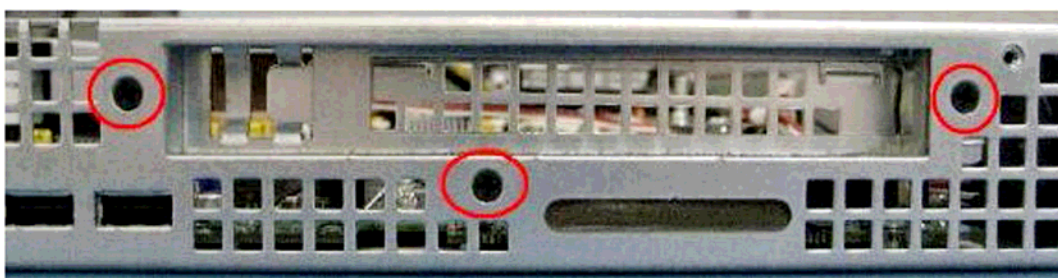


9. Verify that the EMI gasket under the PCI mounting bracket remains securely adhered to the system with no visible damage.

10. Install the PCI cover (FAB-COVERPCI0570) to the rear of the chassis using 3 screws (MSC-M3/.5/5MM/B).




11. Remove and discard the 2 motherboard screws indicated below.

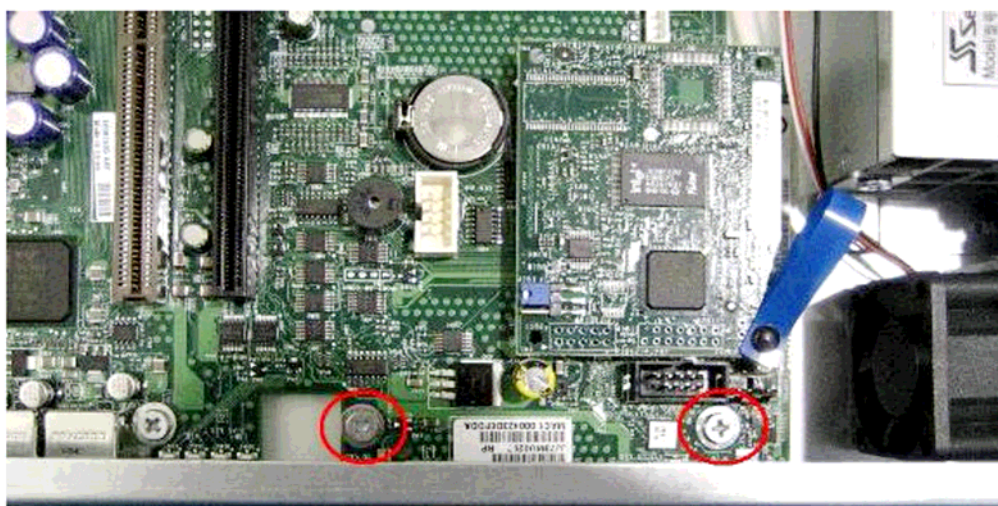


12. Install the 2 standoffs (MSC-STNDFF 3/8-.5) using a 3/8" hex nut driver (MSP-3/8NUTDRIVER) where the 2 motherboard screws were removed.


NOTE: Be careful not to damage any motherboard components while securing the standoffs.

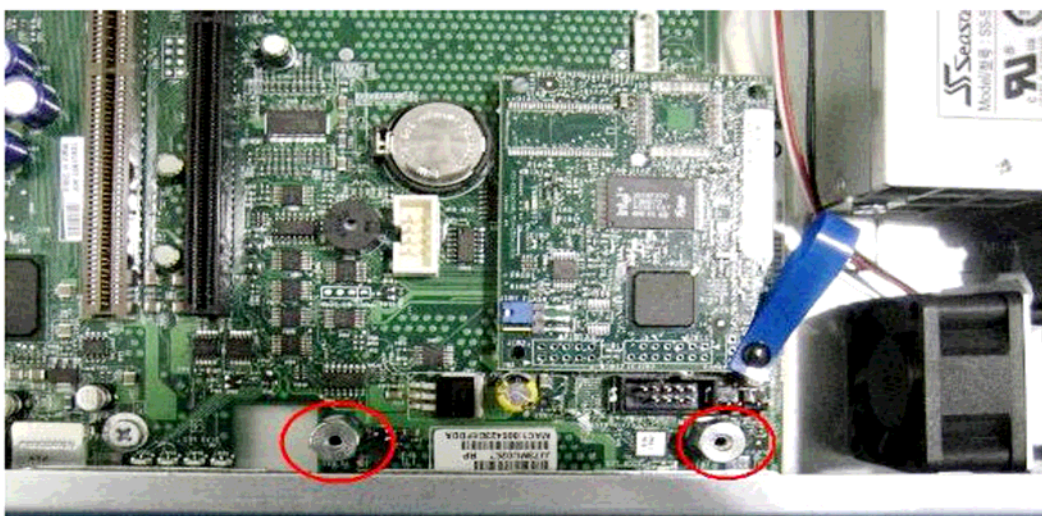

13. Install the VRAC mounting bracket (FAB-BRKT0569) to the standoffs using 2 screws (MSC-6/32/.1875/Z).

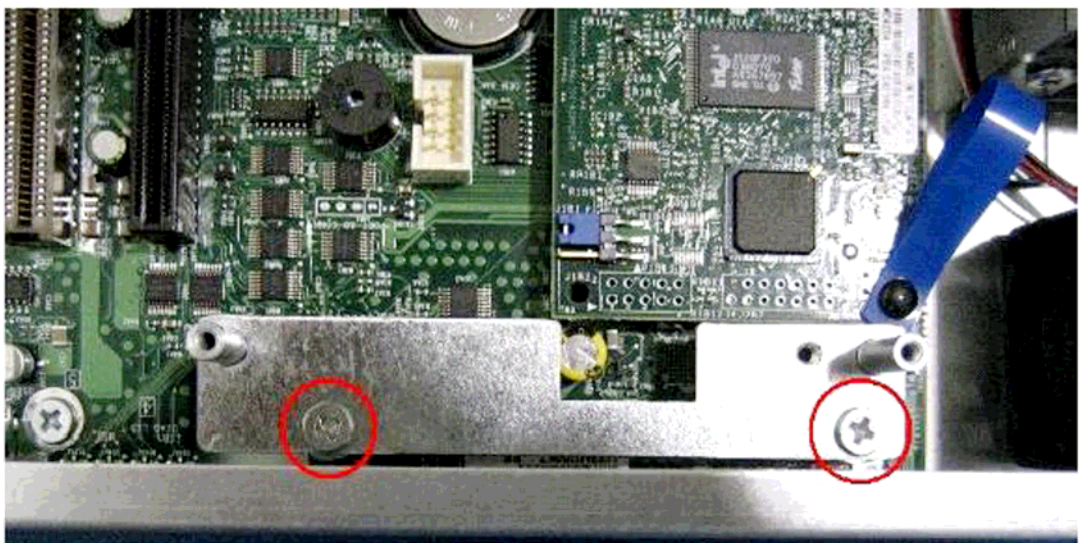


14. Loosen, but do **NOT** remove both screws on the VRAC alignment bracket.

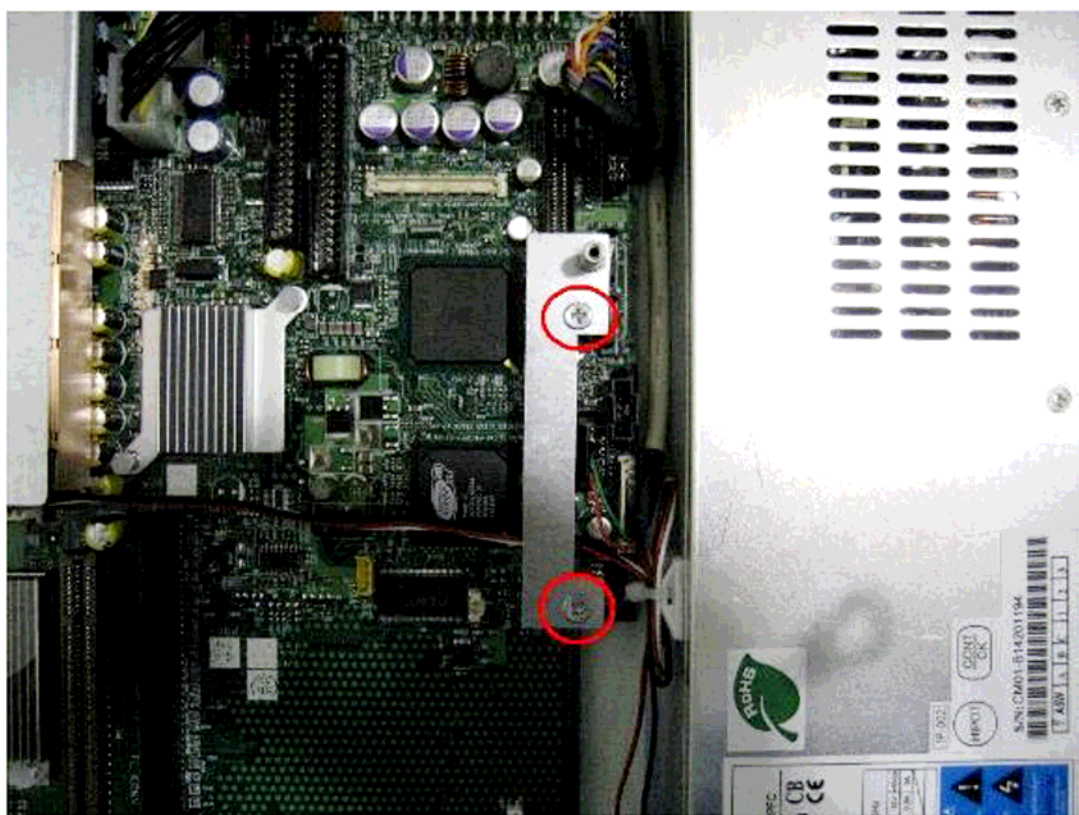



15. Connect the VRAC mounting jig (FXT-JIGVRAC) to the riser that was removed in step 6.

NOTE: To prevent damage to the riser and the motherboard, it is critical to utilize the mounting jig as specified to properly position the alignment bracket!

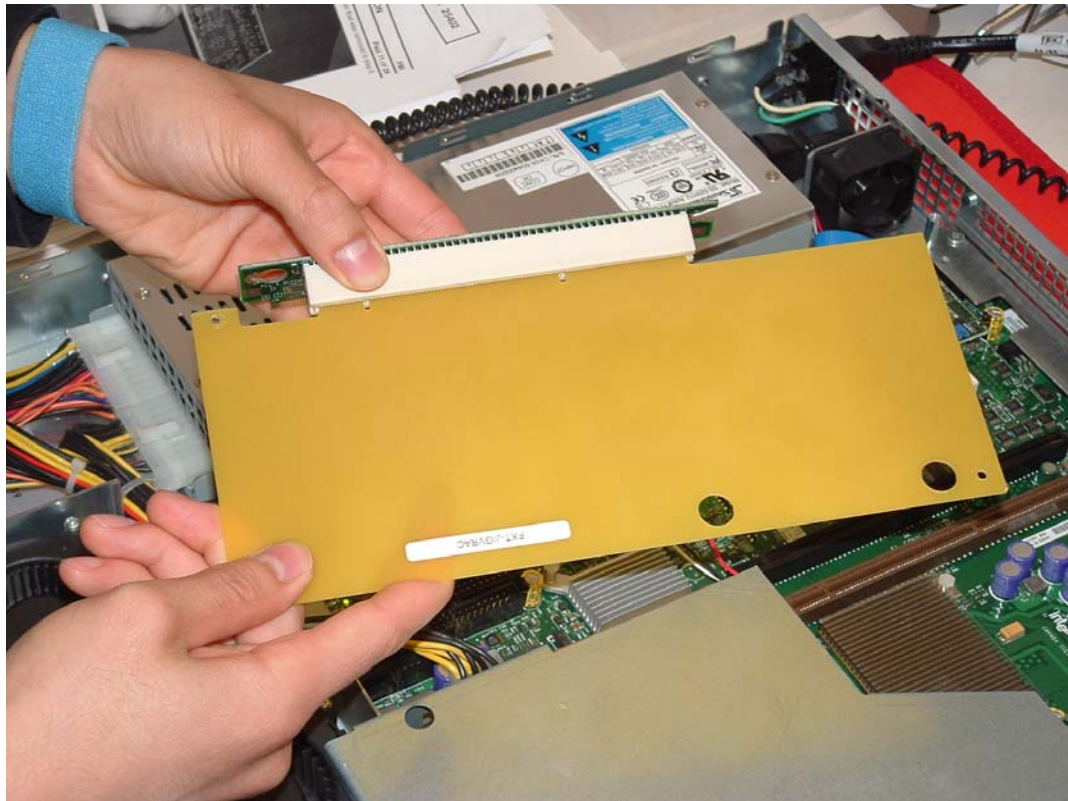


16. Install the jig / riser assembly to the motherboard's black PCIExpress slot by pushing straight down on both ends of the riser.


NOTE: Verify the fan wire is not pinched while installing the jig / riser into the motherboard.

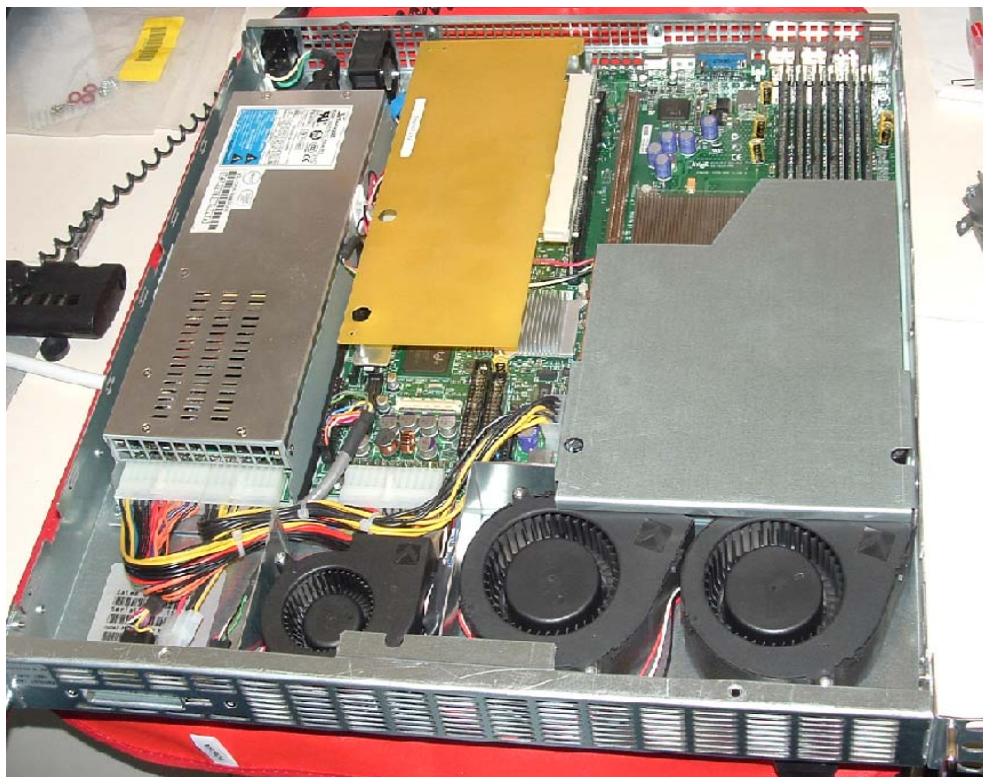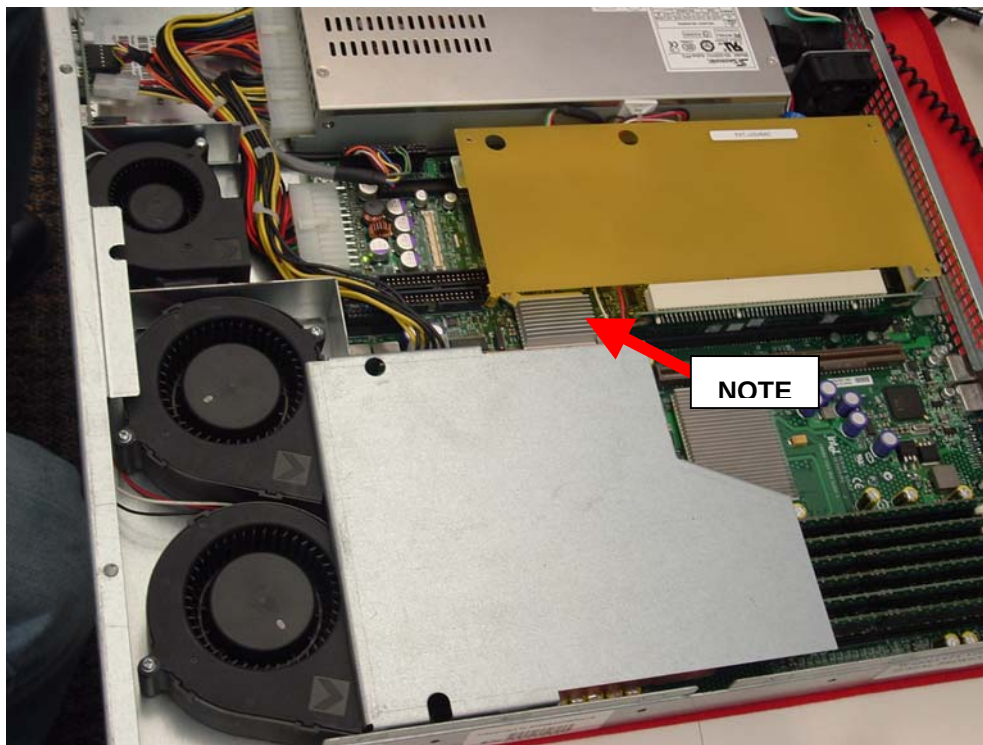


17. Secure the VRAC mounting jig (FXT-JIGVRAC) to the VRAC mounting bracket using 2 screws (MSC-M3/.5/5MM/Z).

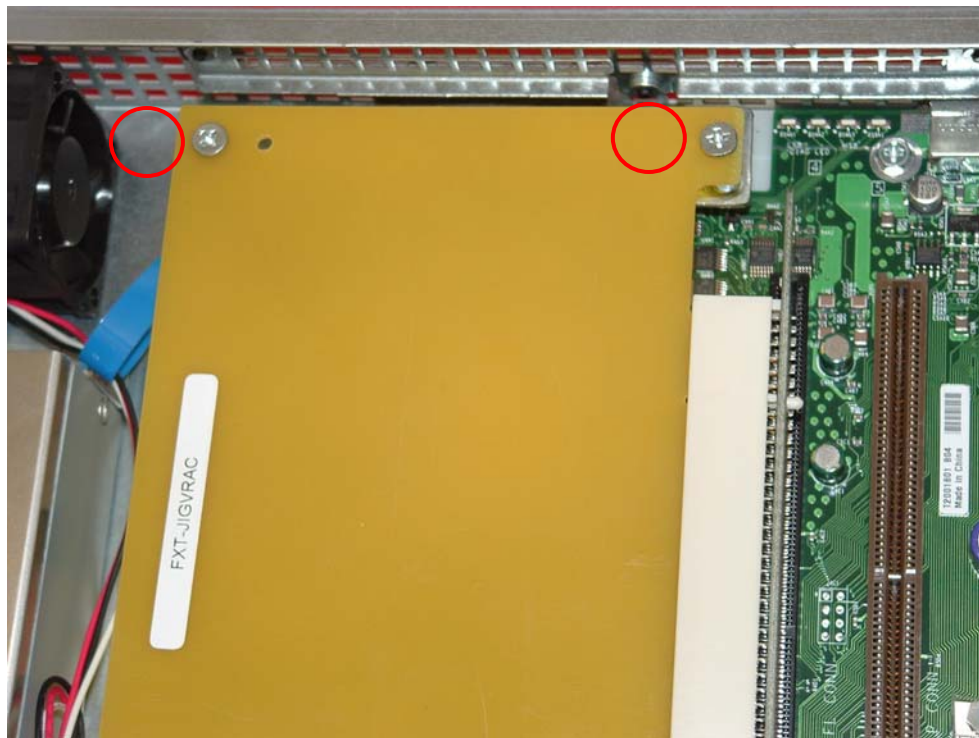


18. Line up the VRAC alignment bracket to the hole on the VRAC card by sliding the bracket.

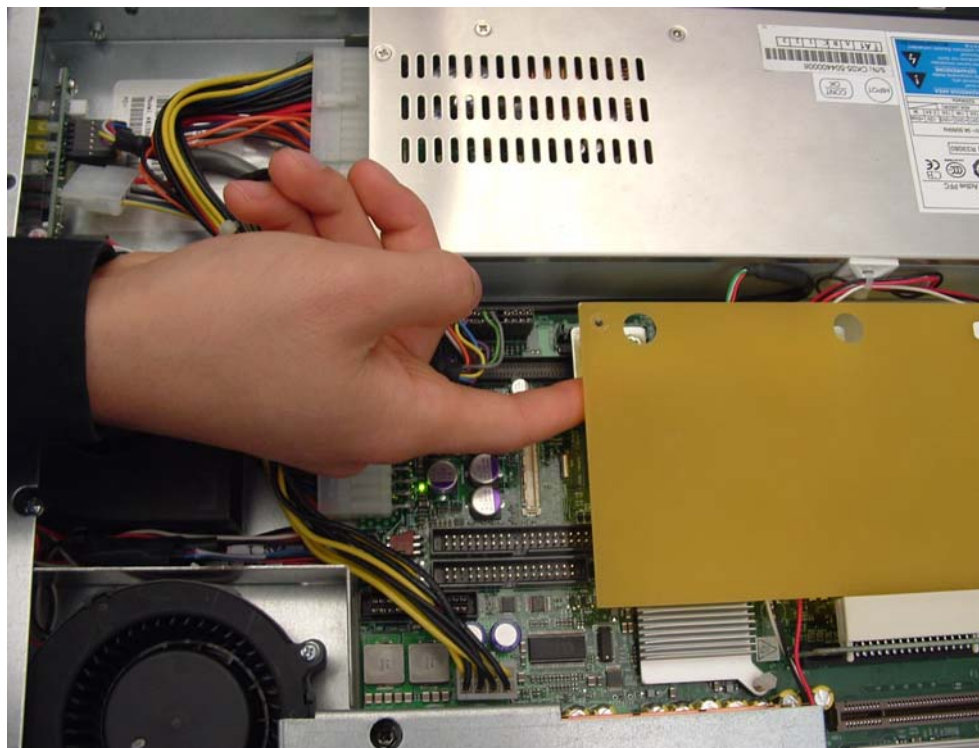


19. Secure the VRAC mounting jig using 1 screw (MSC-M3/.5/5MM/Z).

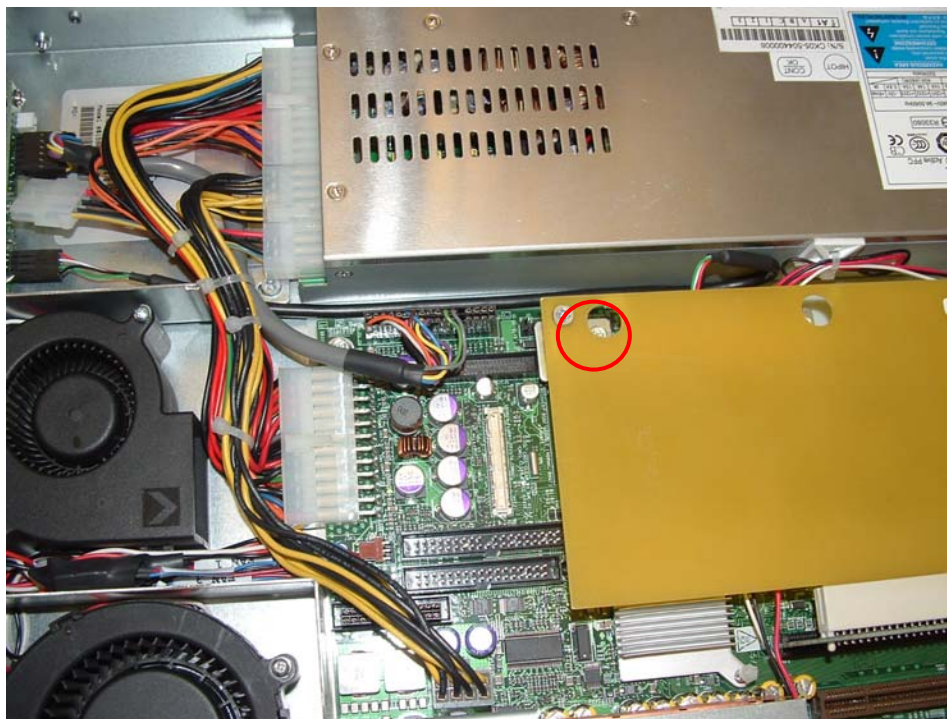

20. Place the driver bit through each hole on the VRAC mounting jig (FXT-JIGVRAC) to secure both screws on the alignment bracket.

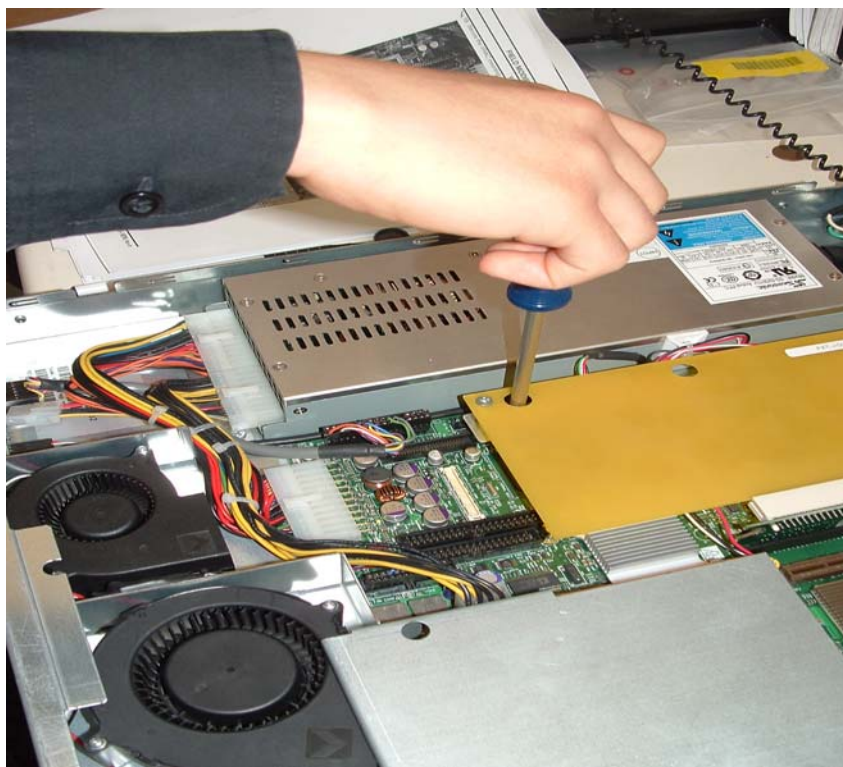


21. Remove all 3 screws that are securing the VRAC mounting jig (FXT-JIGVRAC).

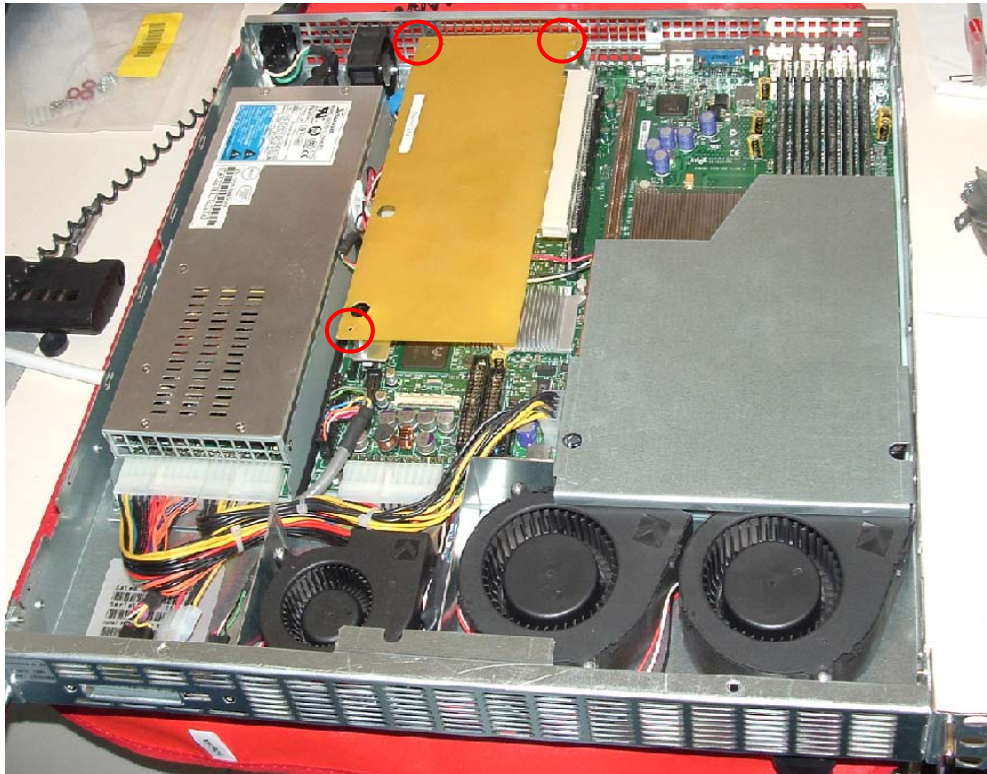


22. Remove the jig / riser assembly from the system by pulling straight up on both ends of the riser.

23. Disconnect the VRAC mounting jig (FXT-JIGVRAC) from the riser.

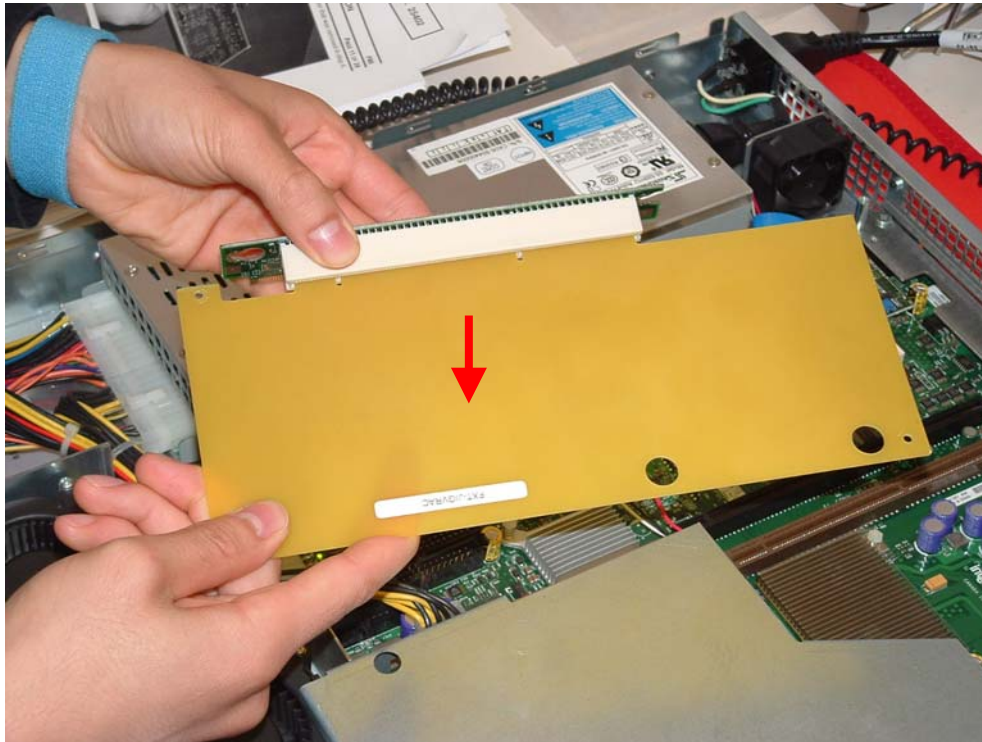


24. Connect the functional VRAC card to the riser.

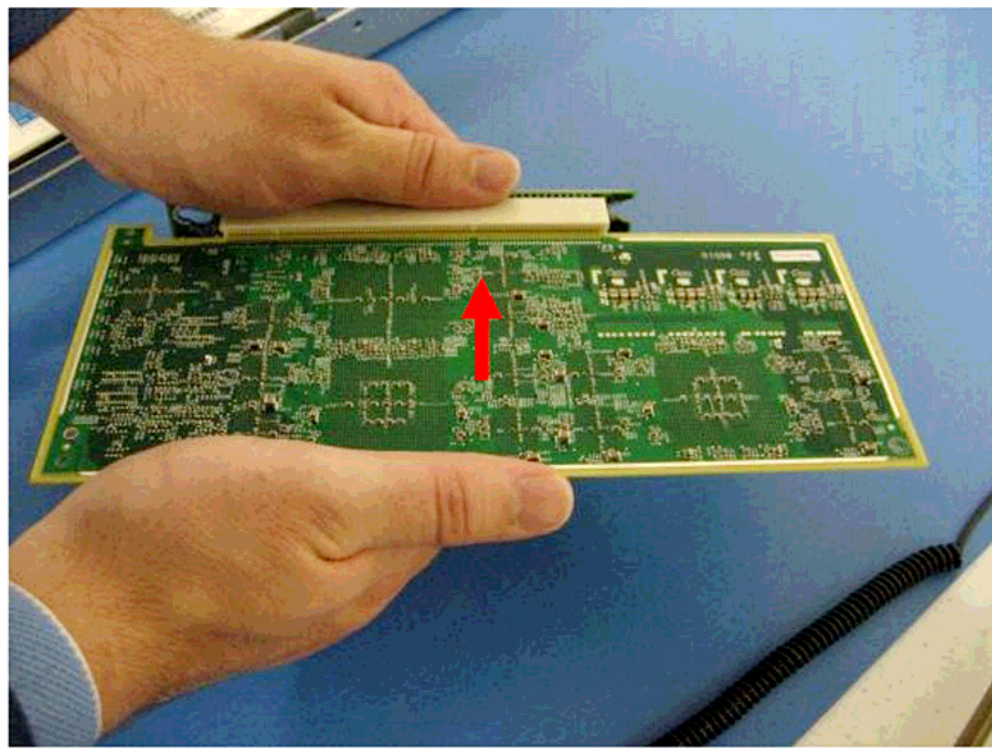


25. Install the VRAC / riser assembly to the motherboard's black PCI-Express slot by pushing straight down on both ends of the riser.

NOTE: Verify the fan wire is not pinched while installing the jig / riser into the motherboard.

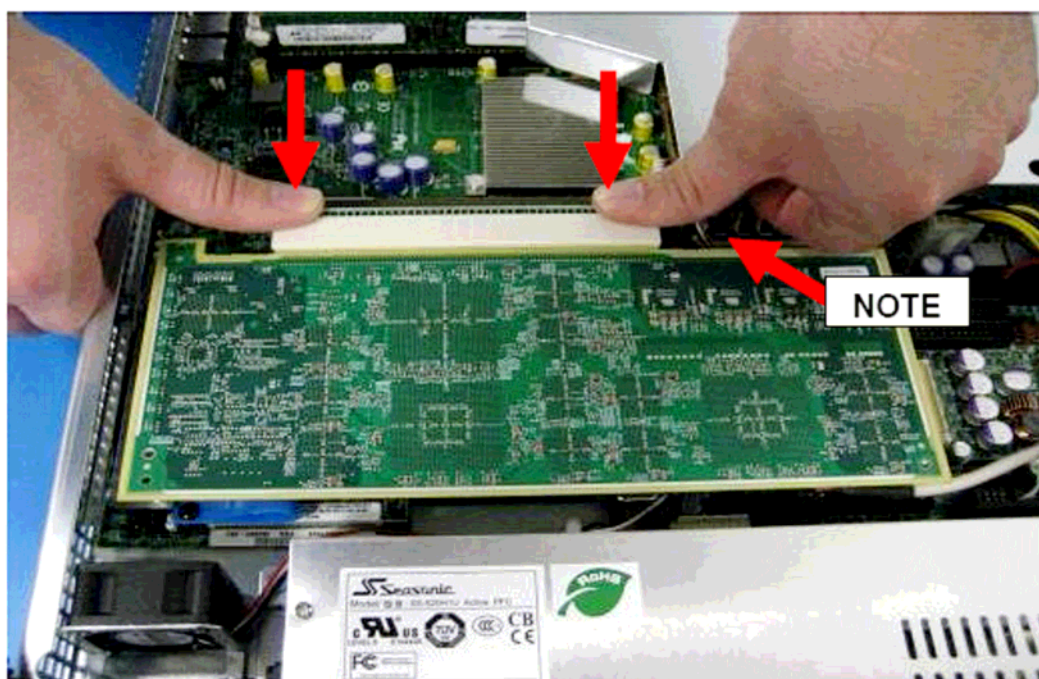


26. Place 1 spacer (MSP-SPACER03125) over the hole on the VRAC card closest to the riser.

CAUTION: Be careful not to drop the spacer inside of the system. Do **NOT** place any other spacers on the VRAC card at this time.

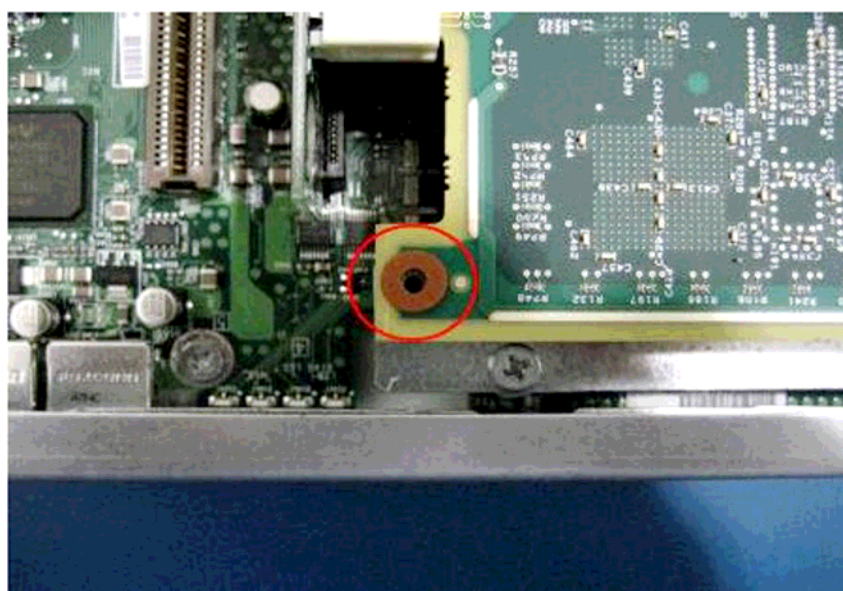


27. Secure the VRAC card using 1 screw (MSC-M3/.5/5MM/Z).

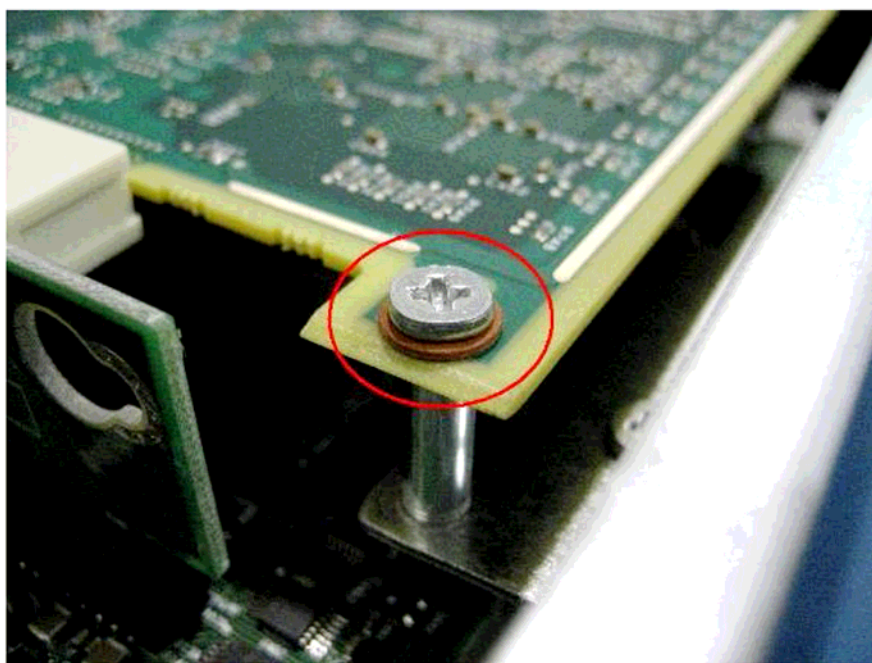


28. Place 1 spacer (MSP-SPACER03125) over the hole on the VRAC card closest to the rear exhaust fan.

CAUTION: Be careful not to drop the spacer inside of the system. Do **NOT** place any other spacers on the VRAC card at this time.

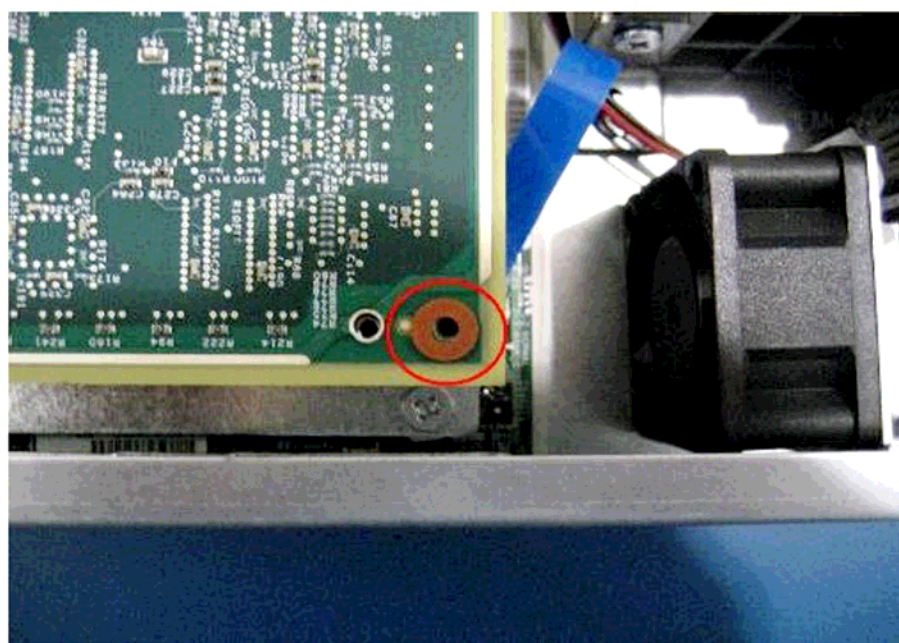


29. Secure using 1 screw (MSC-M3/.5/5MM/Z).

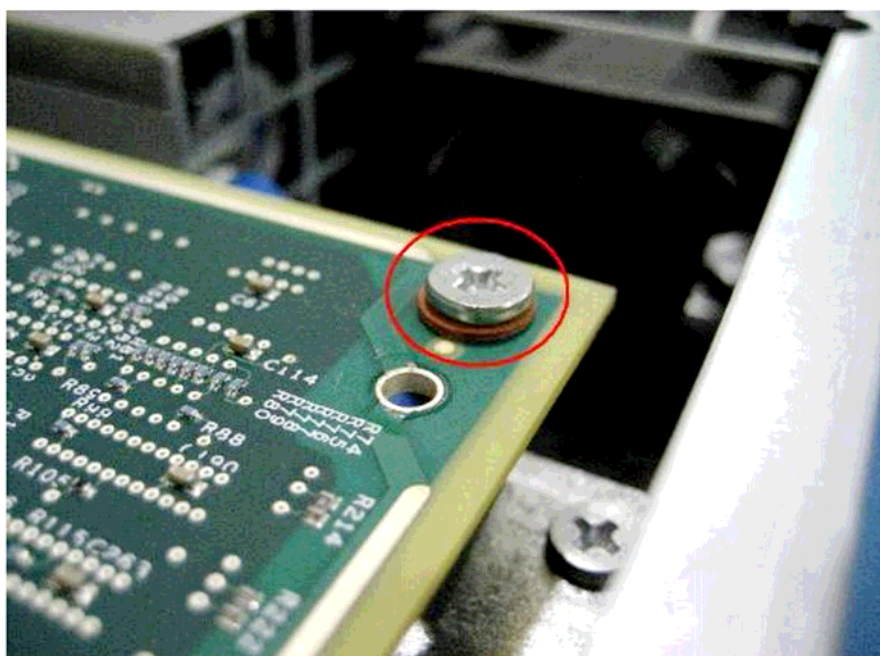


30. Place 1 spacer (MSP-SPACER03125) over the hole on the VRAC card closest to the VRAC alignment bracket.


CAUTION: Be careful not to drop the spacer inside of the system.

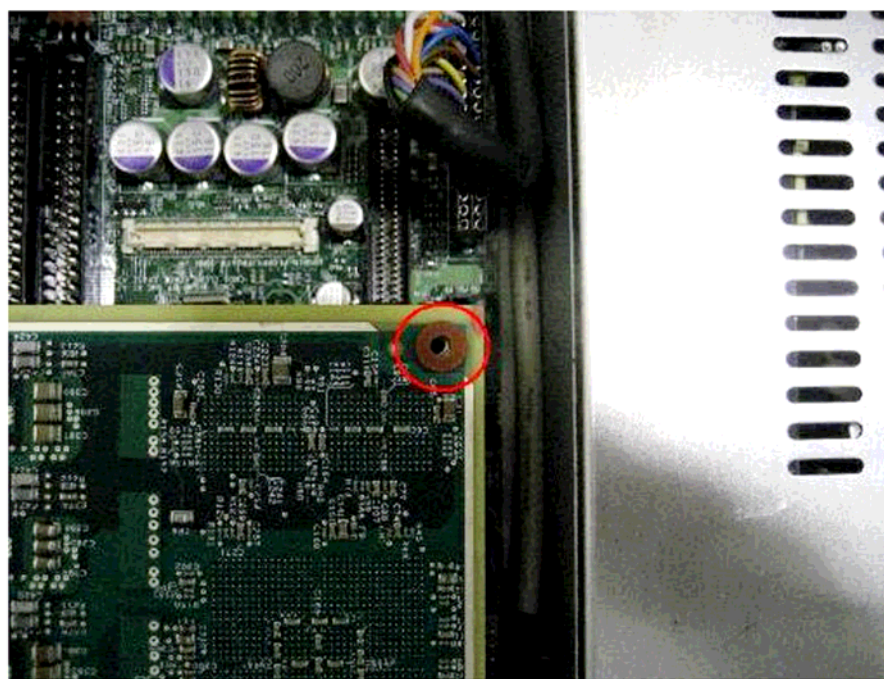

31. Secure using 1 screw (MSC-M3/.5/5MM/Z).

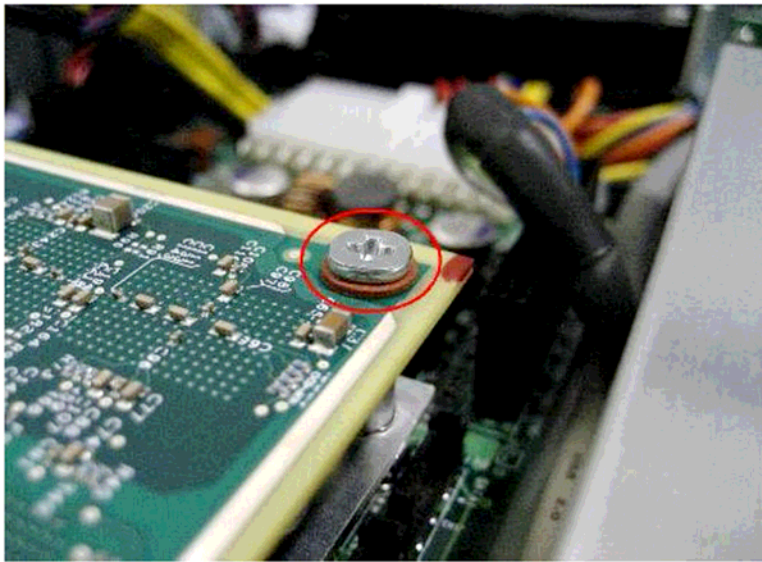


32. Verify that the VRAC card is secured using 3 spacers (MSPSPACER03125) and 3 screws (MSC-M3/.5/5MM/Z).

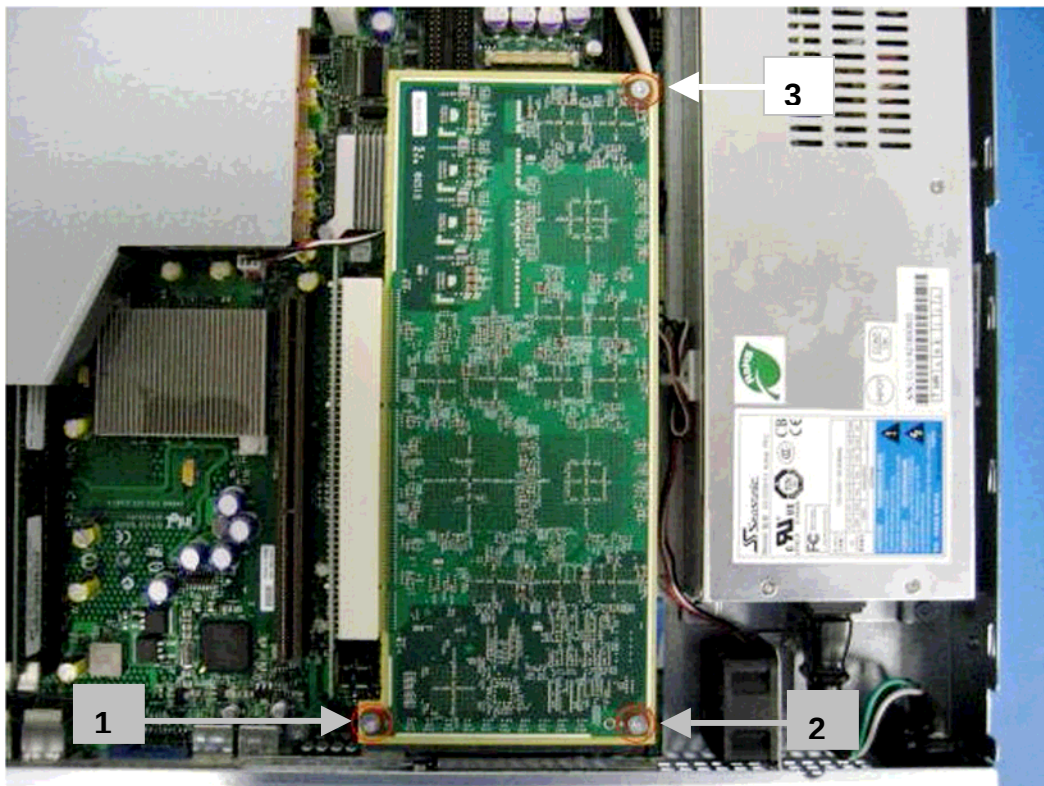

33. Verify that the metal tabs are bent out on both sides of the system as shown. These ensure that the lid makes a good connection with the system.

NOTE: If any of the tabs are not bent out, use a small flat-head screw driver to bend them out.

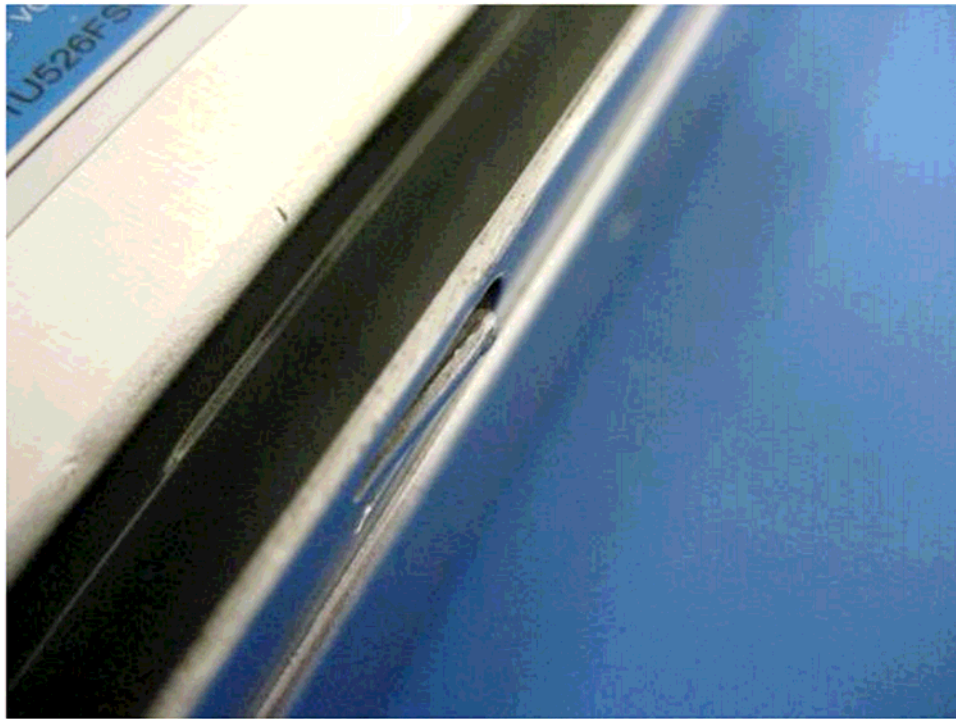

34. Install the lid ensuring that the front tabs are under the front bezel.
35. Secure the lid using 4 Silver Lid Screws (MSC-M3/.5/4/Z).

36. Place the label 46-302200P12 on the side of the system on top of the warranty seal (referenced in Step 1) to indicate that the rework procedure has been completed.




37. Using the Sharpie Marker provided in the kit, locate the label at the back of the VIG. Manual change the model number from 5159834-2 to **5159834-3**. Change VIG-2 (JR) to **VIG-3 (JR)**.




Install the VIG Node into Console

NOTE: If this site has only one VIG Node, the VIG **MUST** be placed in the bottom standard slot 1 position. The Ethernet cable identifying the bottom IG #1 must be correctly connected between the V/DARC Node and the VIG Node. If two Jarrell VIG Nodes are present, the second VIG Node must be placed in slot 2 (middle).

1. Slide the replacement VIG Node(s) into the Operator Console from the front. The bottom slot is the standard slot reserved for VIG Node #1. VIG Node #2 is positioned directly above VIG Node #1. VIG Node #3 is positioned directly above VIG Node #2.
2. Replace the two (2) front mounting screws on the front of the VIG Node(s).
3. Verify the power switch (if present) located at the rear of the VIG Node(s) is set to the ON position.
4. Replace the rear Ethernet cable(s) between the V/DARC and the VIG Node(s).

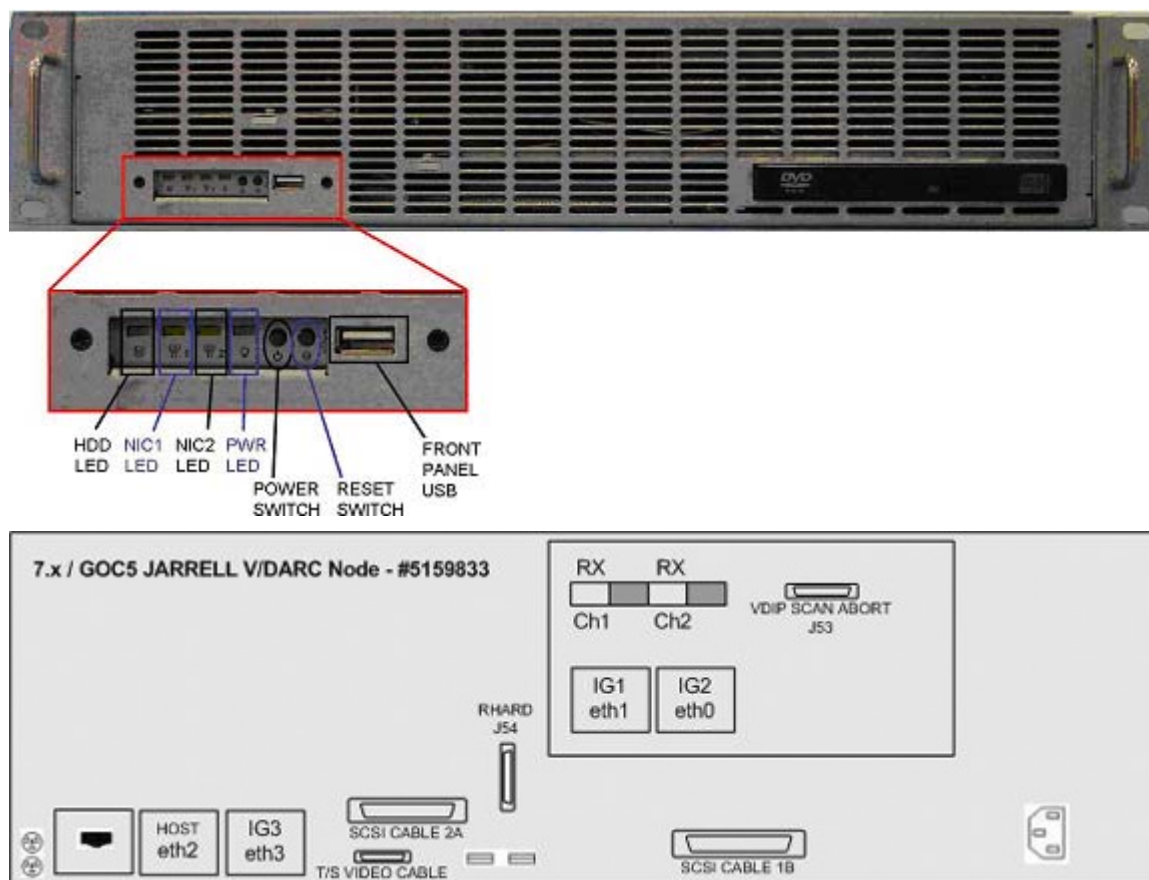

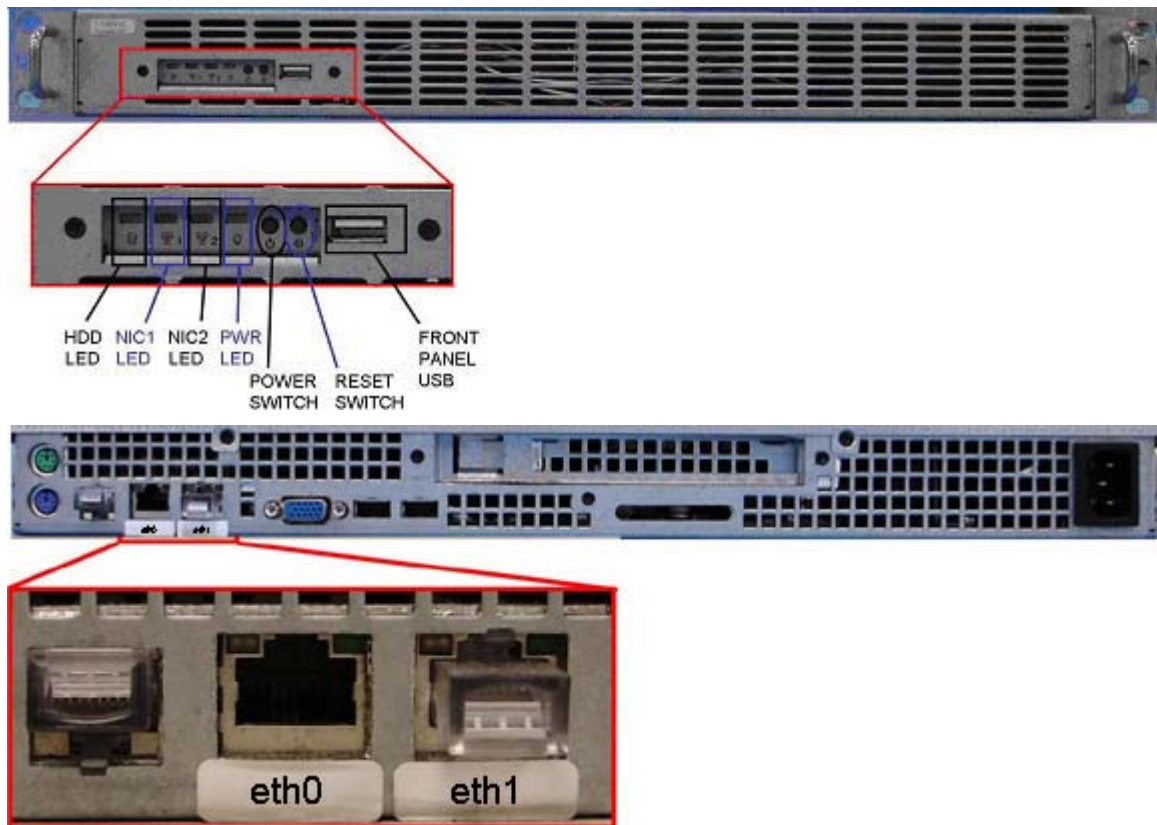

NOTE: This procedure must be performed on all VIGs in the Console. If there are additional VIGs, return to Page 4, Remove VIG Node from Console and repeat the steps for each additional VIG!

Power-ON the Operator Console

1. Power ON the Operator Console at the front panel switch.
2. At the V/DARC Node, verify the rear power switch is set to the ON position.
3. Visually verify the V/DARC Node front panel Power LED is illuminated. If not, press the front panel power switch (3-10 seconds) to apply power to the V/DARC Node.
4. Visually verify the VIG Node front panel Power LED(s) is illuminated. If not, press the front panel power switch (3-10 seconds) to apply power to the VIG Node.
5. At the Disk Array (JBOD) verify two green LEDs at rear panel are ON. If not, press two power switches at rear panel, to apply power to the Disk Array.
6. Stop Applications Software from auto-starting at the 5 second Pink Box.

Verify V/DARC Node Is Up

The V/DARC Node must be up to check the VIG Node(s). Perform `ping` and `rsh` commands on the V/DARC Node, as follows:

1. At the Toolchest, with Application software down, open a Unix Shell.
2. Ping the V/DARC Node:
`{ctuser@hostname} ping darc <ENTER>`
Verify `ping` is successful.
Press `<Ctrl+C>` to stop ping process.
3. Remote Shell to the V/DARC Node:
`{ctuser@hostname} rsh darc <ENTER>`
Verify `ctuser@darc` prompt is present.
4. Exit from the V/DARC Node Remote Shell:
`{ctuser@darc} exit <ENTER>`


NOTE: The Unix Shell will be used again, so do not exit or close at this time.

Verify Replacement VIG Node

In some configurations, the GOC5 Operator Console may have up to three (3) VIG Nodes. In the instructions below, be certain to specify `ig1`, `ig2` or `ig3`, as appropriate for the VIG Node being tested.

NOTE: Application Software must be down.

1. Power OFF the replaced or repositioned VIG Node at the front panel switch.

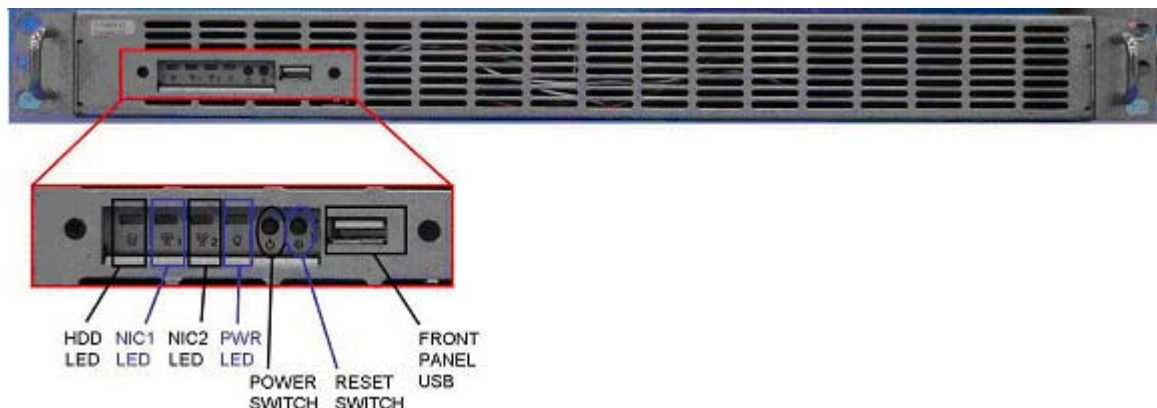

2. Power ON the replaced or repositioned VIG Node at the front panel switch (3-10 seconds).

3. Wait two to four minutes for the VIG Node boot up.

NOTE: Whenever a VIG Node is swapped (physical location or on end of the ethernet cable), the affected VIG Node requires a manual power cycle at the front panel. This manual power cycle is required for the VIG Node to get an IP address.

4. In the Unix Shell, verify the `cat /usr/g/config/INFO` file for number of VIGs active in the Operator Console:

NOTE: If the number of VIG Nodes present is incorrect, then perform the reconfig process as root on the Host with application Software down

```
{ctuser@hostname} cat usr/g/config/INFO <ENTER>
```

5. In the Unix Shell, verify the following IG Node(s) as applicable:

NOTE: For communication problems, refer to ethtool T/S Note in Appendix A.

```
{ctuser@hostname} rsh darc <ENTER>
```

```
{ctuser@darc} exit <ENTER>
```

```
{ctuser@hostname} vrac_flash_update <ENTER>
```

```
answer: y <ENTER>
```

NOTE: If `vrac_flash_update` fails, try the process again. If it still fails verify the V/DARC and VIG Nodes can successfully perform the remote shell commands.

```
{ctuser@hostname} rsh darc <ENTER>
```

```
{ctuser@darc} su - <ENTER>
```

```
Password: #bigguy <ENTER>
```

```
[root@darc] rsh ig1 <ENTER>
```

```
[root@ig1] vrac_menu -a <ENTER>
```

- If there is a second IG, verify diagnostics pass.

```
[root@ig1] exit <ENTER>
```

```
[root@darc] rsh ig2 <ENTER>
```

```
[root@ig2] vrac_menu -a <ENTER>
```


- If there is a third IG, verify diagnostics pass.

```
[root@ig2] exit <ENTER>
```

```
[root@darc] rsh ig3 <ENTER>
```

```
[root@ig3] vrac_menu -a <ENTER>
```

- Exit out of IG / DARC.

```
[root@ig3] exit <ENTER>
```

```
[root@darc] exit <ENTER>
```

```
{ctuser@darc} exit <ENTER>
```

```
{ctuser@hostname}
```

6. Start Application Software

```
{ctuser@hostname} st <ENTER>
```

Verify Recon Data Path Diagnostic

1. Select the **Service** Icon.
2. Select **[DIAGNOSTICS]**.
3. Select **[Recon Data Path]**.
4. Select **[1]** and **[ALL]**
5. Select **[RUN]**
6. Verify Recon Data Path tests pass.
7. Select **[Dismiss]**.

Finishing Up

1. Perform sanity scans (e.g., scout, axial and helical).
2. Verify the image(s) reconstruct and are properly displayed.
3. Install the front and rear Operator Console covers.
4. Remove FMI Instructions and all discarded parts from site.
5. Close dispatch per normal region process.

Appendix A – ethtool T/S Note

For an Ethernet cabling connection issue, use the command '**ethtool -p eth#**' (where # represents the Ethernet address) as root on the VDARC Node or VIG Node to determine if the Ethernet cable is correctly connected from the VDARC Node to the VIG Node. The port will blink or flash when activated and <Ctrl+C> will stop the test. This command can be used on the Host, VDARC Node and the VIG Node(s).

Example for VDARC (Jarrell and Westville) Ethernet cable:

```
[root@dar] ethtool -p eth2 (<Ctrl+C> to stop flashing Jarrell Ethernet port LED)
```

```
[root@dar] ethtool -p eth4 (<Ctrl+C> to stop flashing Westville Ethernet port LED)
```

Example for VIG Node Ethernet cable:

```
[root@ig1] ethtool -p eth0 (<Ctrl+C> to stop the flashing Ethernet port LED)
```

```
[root@ig2] ethtool -p eth0 (<Ctrl+C> to stop the flashing Ethernet port LED)
```


```
[root@ig3] ethtool -p eth0 (<Ctrl+C> to stop the flashing Ethernet port LED)
```

Appendix B – Verifying Jarrell VIG

This section is to help verify that the VIG on site is a Jarrell VIG as opposed to a Westville VIG.

Hardware Verification

The best way to tell the difference between a Jarrell and Westville VIG is to examine the front panel power switch. Below is a picture of the front panel of 3 Jarrell VIGs in a console. Note, the power switch is red, with a series of horizontal LED lights.

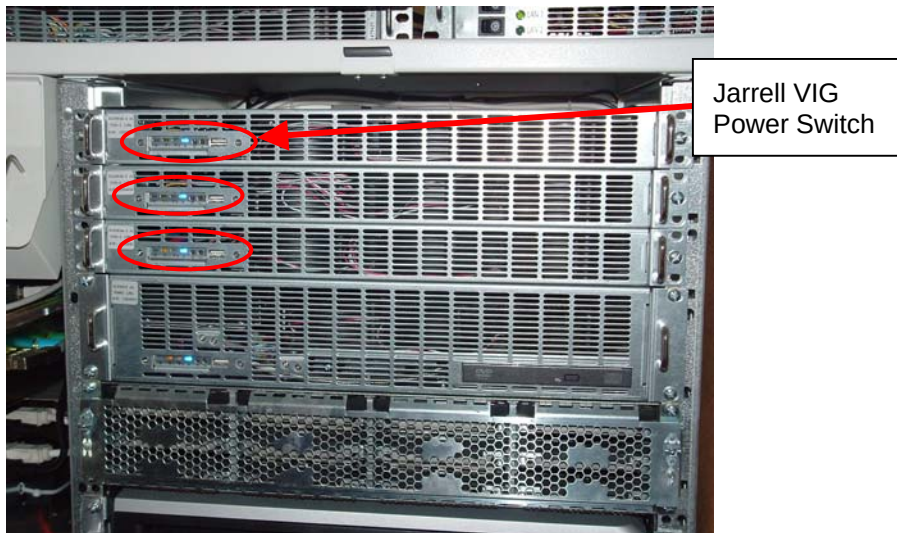

Below is a picture of a Westville VIG for reference. Note, the power switch is a black toggle switch, with LEDs aligned vertically to the right. This FMI does **NOT** apply to Westville VIGs.



Remote Verification

Verification of the type of IG(s) installed can be done on site, or remotely using the following commands:

IN A UNIX SHELL, OR REMOTELY.

1. Type: **su -** <ENTER> (*become root*)
2. Enter root password. <ENTER>
3. Type: **rsh darc** <ENTER> (*rsh to darc, then igs*)
4. Type: **rsh ig1** <ENTER>

(Output will look similar to below)

connect to address 10.0.1.2 port 543: Connection refused

Trying krb4 rlogin...

connect to address 10.0.1.2 port 543: Connection refused

trying normal rlogin (/usr/bin/rlogin)

5. Type: **dmidecode | less** <ENTER>

(Output will look similar to below)

dmidecode 2.6

SMBIOS 2.3 present.

56 structures occupying 2431 bytes.

Table at 0x000FD4A0.

Handle 0x0000

DMI type 0, 20 bytes.

BIOS Information

Vendor: Intel Corporation

Version: SE7520**JR**23.86B.P.10.00.0087.120820051348

Release Date: 12/08/2005

Address: 0xF0000

Runtime Size: 64 kB

ROM Size: 2048 kB

Characteristics:

PCI is supported

PNP is supported

BIOS is upgradeable

BIOS shadowing is allowed

ESCD support is available

Look at Version Line,
JR means Jarrell
WV means Westville,
this one is a Jarrell!

Westville will look like this:

Version: S**WV**25.86B.0218.P28.0405111912

Boot from CD is supported
 Selectable boot is supported
 EDD is supported
 3.5"/720 KB floppy services are supported (int 13h)
 3.5"/2.88 MB floppy services are supported (int 13h)
 Print screen service is supported (int 5h)
 8042 keyboard services are supported (int 9h)
 Serial services are supported (int 14h)
 Printer services are supported (int 17h)

6. Type: **q <ENTER>** (type q to quit at the colon)
7. Type: **logout <ENTER>** (returns to darc)
 - a. At this point, rsh to IG2 and IG3 as needed, or continue with step 8.
8. Type: **logout <ENTER>** (exits darc, returns to host)
 rlogin: connection closed.

NOTE: When rsh to ig, if you get the following: "No route to host" then either IG is down or there is no IG. Remember, there can be 1, 2, or 3 IGs.

```

[root@darc ~]# rsh ig2
connect to address 10.0.2.2 port 543: No route to host
Trying krb4 rlogin...
connect to address 10.0.2.2 port 543: No route to host
trying normal rlogin (/usr/bin/rlogin)
ig2: No route to host
  
```

9. To verify the number of IGs in a system, use the following command:
 - a. Type: **cd /usr/g/Config <ENTER>**
 - b. Type: **more INFO <ENTER>**

```

setenv SERVICE_ID 123456VCT
setenv HOSPITAL_NAME "xxxxxxx Medical Center"
setenv MENUTZ "America"
setenv SUNTZ "US/Central"
setenv TZ "CST6CDT US"
setenv TZTEXT "Chicago"
setenv REGENERATE_DATABASE no
setenv NEXT_MOD_NUMBER 1
          
```

FMI	FIELD MODIFICATION INSTRUCTION	FMI
25402	VCT VIG EMC Update	PAGE 34 OF 35
	setenv NEXT_DIAGNOSTIC_EXAM_NUMBER 50000 setenv GATEWAY_HOSTNAME ct01 setenv GATEWAY_IP 10.46.65.58 setenv GATEWAY_NETMASK 255.255.255.0 setenv GATEWAY_BROADCAST 10.46.65.0 setenv GATEWAY_DEFAULT 10.46.65.1 setenv DEV_ENV no setenv DEVELOPMENT_SUBNET 192.9.220 setenv TUBETYPE 16 setenv PDUTYPE NGPDU setenv DASMTYPE 2 setenv KV_VALUES "80 100 120 140" setenv DEFAULTARCHIVE DICOM setenv ADVANTAGERW FALSE setenv CTIDATETIMEFORMAT 4 setenv TIMESEPAR / setenv WEIGHT lb setenv DOCTITLE "radiologist" setenv CAMERATYPE "" setenv SUITEID CT01 setenv FILTERTYPE 2 setenv CAMERA no setenv DEFFORMAT 12 setenv DOGREYSCALE yes setenv DOSLIDES no setenv DOZOOM no setenv FILMINTERP sharp setenv PRINTER_HOSTNAME setenv PRINTER no setenv O2_PRESENT n setenv DASTYPE VDAS64 setenv EKG_CONSIPADDR 10.44.22.20 setenv AWDCON no setenv AWDCIPADDR 10.44.22.22 setenv INROOMSCAN off setenv GANTRY_TYPE vct_40 setenv RECONHARDWARE GRE	

setenv REV_NO 1
setenv GANTRYLAYOUT GantryOnLeft
setenv LANGUAGE en_US
setenv AV_LANG 1
setenv DOSE y
setenv TABLETYPE GT_LONG

setenv NUMIG 3

See line that indicates NUMIG,
in this case, there are 3.

setenv IGNUMLIST 1|3
setenv CONSOLETYPE GOC5
setenv CONSOLELIST GOC5
setenv GANTRYSUBTYPE LS
setenv HRCTTITLE ct01
setenv HRPORTNUM 1003
setenv HRAETITLE NovaWorklist
setenv HRIPADDR 10.46.65.100
setenv IIP_SERVICE_ID 123456VCT