

**ADMINISTRATIVE INFORMATION****FMI No: 25405****FMI Title:** Service Pack 1.5 for CT Software Releases 07MW11.10, 07MW18.4, 07BW08.3 and 07BW08.4 for non-English Speaking Sites

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS – AM:	Issue Date:
☐ Mandatory	GEMS – E:	
☐ Mandatory on Request	GEMS – A:	Close Date:

- ☒ **Staged FMI** **Customer List** provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
- ☒ **Parts** and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
- ☐ **Parts Ordering:** GEMS - Americas, Contact your Region logistics coordinator
GEMS - Europe, Order FMI kits from GPO/EPL Buc
GEMS - Asia, Contact your region logistics coordinator
- ☐ **No Parts** are required to complete this FMI.
- ☐ **Batch with** _____
- To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
- ☐ **Non-Staged** **No Customer List is provided.** Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS – Americas	GEMS - Europe	GEMS – Asia
GE Medical Systems On-Line Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:** Multiple – Refer to “Effectivity” section**FMI Instructions Part No.:** 5272905-1EN, Rev 1**Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS – Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25405

FMI Title: Service Pack 1.5 for CT Software Releases 07MW11.10,
07MW18.4, 07BW08.3 and 07BW08.4 for non-English
Speaking Sites

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments _____

Effectivity

This FMI is deployed to a limited group of CT-only systems made up of LightSpeed 7X, 5X, RT16, RT4, BrightSpeed 4/8/16-slice, and BrightSpeed Select products. This “limited group” is further defined as non-English speaking sites that were part of the following product family shipments:

- LightSpeed 7X systems, FMI 25394’s F3-pilot and early forward production units where 07MW18.4 was cut-in, but not Service Pack 1.5 software (manufacture date range 8-Jun-07 till 23-Aug-07)
- LightSpeed 5X systems manufactured with one of the software releases listed in the “Prerequisites” section of this FMI instruction, but not Service Pack 1.5 software (manufacture date range 1-Aug-07 till 23-Aug-07)
- RT16 & RT4 systems manufactured with one of the software releases listed in the “Prerequisites” section of this FMI instruction, but not Service Pack 1.5 software (manufacture date range 8-Jun-07 till 26-Sep-07)
- BrightSpeed & BrightSpeed Select systems manufactured with one of the software releases listed in the “Prerequisites” section of this FMI instruction, but not Service Pack 1.5 software (BSD Select manufacture date range 20-Jun-07 till 21-Aug-07)

NOTE: DO NOT INSTALL this FMI on any PET/CT systems.

Purpose

This FMI deploys Service Pack 1.5 (SP1.5) software to fix a language translation issue with VoxTool / Reformat screens seen by the operator.

With SP1.5 installed, VoxTool / Reformat screens are properly translated into the language for which the system has been configured.

If SP1.5 software is not installed, the operator will see a mix of English and the selected language for which the system has been configured in VoxTool / Reformat user interface screens.

Prerequisites

The CT-only system must be loaded with one of the software releases listed below, and its associated service pack software as listed, before installing this FMI's Service Pack 1.5 software:

- 07MW11.10, deployed by FMI 25389,
- 07MW18.4, deployed by FMI 25394,
- 07BW08.3 with Service Pack B (p/n 5220681) for RT16, deployed by FMI 22915, or
- 07BW08.4.

Furnished Materials

The following items are included in FMI 25405, kit part number 5272905, Rev 1:

Item	Name	Rev	Qty	Description
1	5272905-1EN	1	1	FMI 25405 Instruction – Service Pack 1.5 for 07MW11.10 & 18.4 and 07BW8.3 & 8.4
2	5234959-15	3	1	Service Pack 1.5 CDROM for software releases 07MW11.10, 07MW18.4, 07BW8.3, and 07BW8.4 with Product Locator Card (PLC)

Tools Required

None

Personnel Requirements

Travel: 1 Hour


Labor: 1 Hour

Field Engineer: 1

Product Safety Fixes

- Mixed-language text in VoxTool / Reformat user interface screens (English and the selected language per the system's configuration)
 - o iTrak PSR 13124499: German Reformat UI screens contains English
 - o iTrak PSR 13125944: French Reformat UI screens contains English

Process Flow

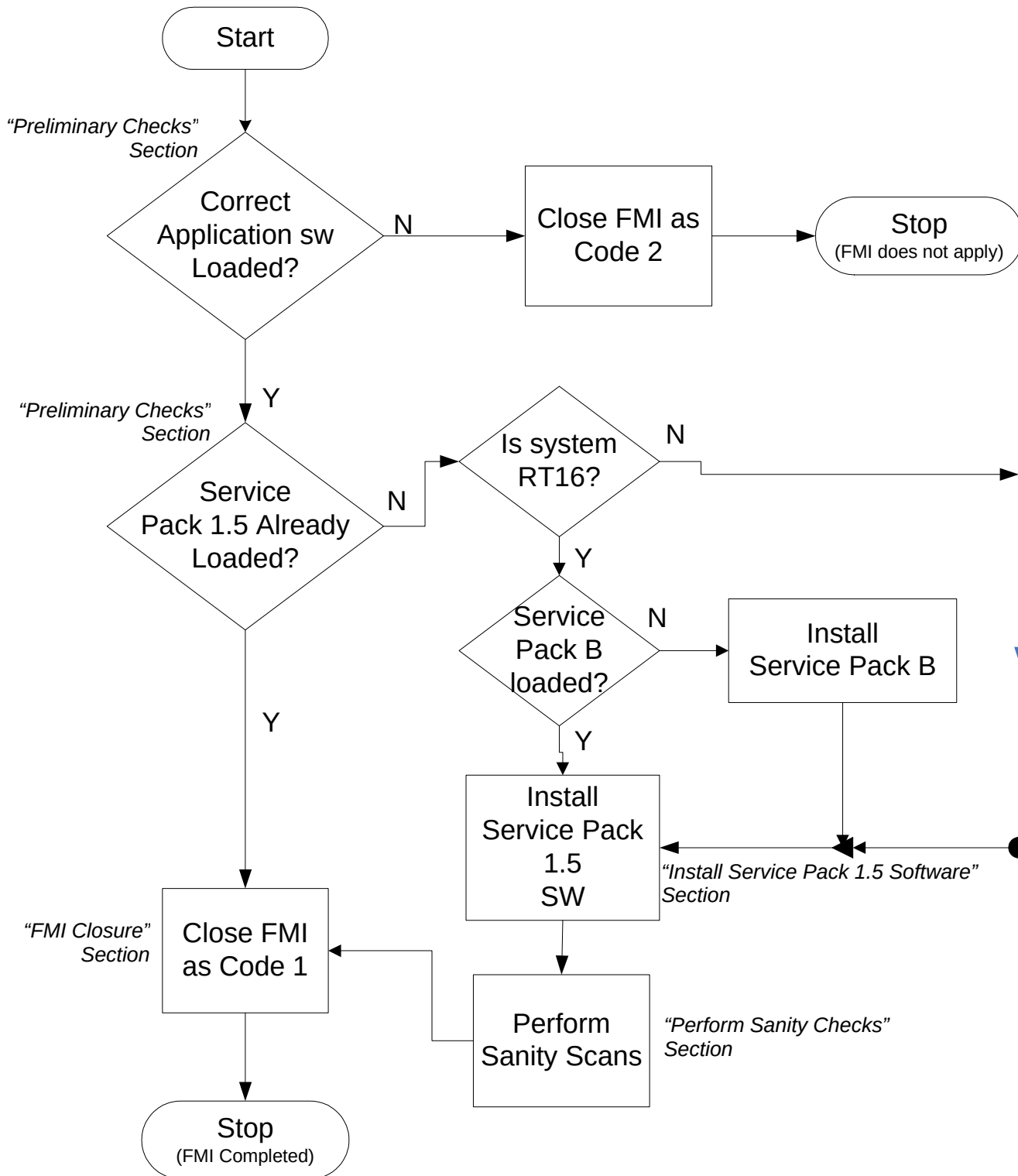

Procedure

Preliminary Checks

Check the software release currently loaded on the host console:

1. Enter Service/Common Service Desktop.
Click "Service" icon on right monitor of host console.
2. View "Home" tab, "Current System Status" area.
3. View "Application Software" field, "Information" column.

Look at the beginning of the software release string for one of the software releases listed in the "Prerequisites" section of this instruction.

STOP if your system does not have one of the software releases listed in the "Prerequisites" section already loaded → this FMI is not to be installed on your system. Close this FMI dispatch as Code 2 and return the FMI kit.

CONTINUE if your system is loaded with one of the software releases listed in the "Prerequisites" section of the FMI instruction → this FMI is to be installed on your system.

4. Check if Service Pack 1.5 software is already loaded.
 - a) Open Unix shell
 - b) Login: **ctuser** [ENTER]
 - c) Type: **showprods | grep -i servicepack** [ENTER]

Example of returned string:

```
07MW11.10_07MW18.4_07BW08.4_GKHCServicePack1.5 08/20/2007  
07MW11.10/07MW18.4/07BW08.4 Service Pack 1.5 ver 1 rel 1.2
```

5. If the returned string does include "Service Pack 1.5," then jump to the section entitled "FMI Closure" because the service pack software is already installed.

If the returned string does not include "Service Pack 1.5," then continue to the section entitled "Manage Earlier-Installed Service Pack Software."

Manage Earlier-Installed Service Pack Software

For all systems, "Service Pack 1" must be removed before installing "Service Pack 1.5." For RT16 systems only, "Service Pack B" software must be installed before installing "Service Pack 1.5."

6. Check if Service Pack 1 is installed.

Review the results of the **showprods** command executed in Step 4c.

- a) If the results show “07MW11.10_servicepack 1 ver 1 rel 1.1” is installed, then remove Service Pack 1 before installing Service Pack 1.5, else its installation will fail.

Remove Service Pack 1 as follows:

- i. Open Unix shell.
 - ii. Type: **cleanMon** [Enter] to shutdown applications.
 - iii. Login as root.
 - iv. Type: **patch_install -u** [Enter] to uninstall patch.
 - v. Select the patch to uninstall from the system by typing **Y** (yes) when queried for which patch to remove.
 - vi. Type: **reboot** [ENTER]
- b) If the results show “07MW11.10_servicepack 1 ver 1 rel 1.1” is not installed, then continue with the next step below.

7. For RT16 systems only, check if “Service Pack B” for 07BW08.3 is installed, p/n 5220681.

- a) Type: **showprods | grep -i 07BW08_SPB** [ENTER]
- b) If the results show 07BW08_SPB is installed, then continue to the next section entitled “Install Service Pack 1.5 Software.”
- c) If the results show 07BW08_SPB is not installed, then load “Service Pack B” now.

Service Pack B shipped inside FMI 22915.

Install Service Pack B as follows:

- i. Open Unix shell.
- ii. Type: **cleanMon** [Enter] to shutdown applications.
- iii. Login as root.
- iv. Place the Service Pack B CD, p/n 5220681, in the host computer’s DVD drive.
- v. Type: **patch_install -c** [Enter] to install “07BW08_SPB-1-1.i386.rpm.”
- vi. Type: **reboot** [ENTER]

Install Service Pack 1.5 Software

Service Pack 1.5 software must be installed for VoxTool / Reformat user interface screens to be properly translated into the language for which the system has been configured.

8. Install Service Pack 1.5 software, CDROM p/n 5234957-15.
 - a) Retrieve Service Pack 1.5 CDROM from this FMI kit.
(Inside software jewel case p/n 5234959-15)
 - b) Insert the CDROM into the HP PC DVD drive.
 - c) Shutdown Apps → open Unix shell and type **cleanMon** [ENTER]
 - d) In Unix shell, Login: **su** – [ENTER]
 - e) Password: **#bigguy** [ENTER]
 - f) Type: **patch_install -c** [ENTER]
 - g) Type: **sync** [ENTER]
 - h) Type: **sync** [ENTER]
 - i) Type: **reboot** [ENTER]
 - j) Remove the Service Pack 1.5 CDROM from the drive and store on-site.
9. Confirm Service Pack 1.5 has been installed:
 - a) Open Unix shell
 - b) Login: **ctuser** [ENTER]
 - c) Type: **showprods | grep -i servicepack** [ENTER]
 - d) If SP1.5 successfully installed, view “ServicePack1.5” in the returned string.

Example of returned string:

```
07MW11.10_07MW18.4_07BW08.4_GEHCSERVICEPACK1.5 08/20/2007
07MW11.10/07MW18.4/07BW08.4 Service Pack 1.5 ver 1 rel 1.2
```

NOTE: If SP1.5 failed to install, refer to Appendix A entitled “Service Pack 1.5 “patch_install” Fails to Mount Device” for an alternate installation procedure.

10. Ensure the system is configured for the language desired by the site.

If the correct language is not set, reconfigure the system now to select the proper language desired by the site (reconfig → preferences → select language → save/exit).

NOTE: If SmartView is used, do not reconfigure the language to French; leave system configured for English.

If reconfigured to French, an issue will arise, SmartView will not display the images and ImageWorks must be used to view them.

Perform Sanity Checks

Perform sanity checks to confirm that Service Pack software has been installed correctly (no mixed languages on operator screens) and to ensure the system is functioning properly before returning it to the customer for patient imaging.

11. Perform the following sanity scan and image review:

a) Prescribe Abdominal or Head exam

- 50 images
- 0.5 secs rotation
- 120kV, 40mA

b) Execute scan.

c) Select ImageWorks.

d) Select newly-created image series.

e) Select Reformat.

f) Confirm the model is built and images displayed.

g) Confirm operator screens do not show English text mixed with the configured language

(click "Filming Tools" → click "Batch Film" → ensure no mixed text).

The system must present operator screens without mixed language text, scan, and display images without error before returning the system to the customer.

FMI Closure

12. Remove FMI Instructions from the customer site (do not leave them on site).

13. Close FMI per regional procedure and close.

Appendix A

Service Pack 1.5 “patch_install” Fails to Mount Device

Installing service pack software using the “patch_install” command may fail because the “/dev/cdrom” device (for the host or SCSI tower DVD drive) may not successfully mount. This is an intermittent issue associated with operating system OS 6.2.9 (the OS used with CT-only software releases 07MW11.10, 07MW18.4, 07BW08.3, 07BW08.4).

The following procedure may be used to mount the drive and load service pack software.

1. Open a Unix Shell

2. Verify Application software is down.

If Application software is not down, perform the following script.

```
{ctuser@hostname}    cleanMon
```

3. Wait approximately 1-2 minutes for Application software to go down.

4. Become root:

```
Login:      su -
```

```
Password:  #bigguy
```

5. As root on the Host type:

```
[root@hostname ~]#    start_udev
```

6. Wait at least 15 seconds after the OK.

7. Retry the Service Pack installation command.

8. If this fails, retry the **start_udev** command, then Service Pack 1.5 installation, several more times.

Recall, this is an intermittent issue (cdrom mounting correctly). Multiple attempts can yield success.

9. If the retry of the command continues to fail, perform a reboot as root
[root@hostname ~]# **reboot**

NOTE: When reboot is performed as “root,” applications software does not auto-start.

Appendix B

How to Uninstall Service Pack 1.5

The following instructions are to be executed to uninstall Service Pack 1.5 (Service Pack 1.5 CDROM, p/n 5234957-15, is not needed to perform its de-installation):

1. Open Unix Shell.
2. Type: **cleanMon** [ENTER] to Shut down all CT Applications.
3. Login as root.
4. Type: **patch_install -u** [ENTER]
5. Select the patch to uninstall from the system by typing **Y** (yes) when queried for which patch to remove.
6. Type: **reboot** [ENTER]