GE Healthcare

1.5T General Purpose Flex Coil

Service Manual

Direction 2336152
Rev 6
© 2002, General Electric Company,
All Rights Reserved

TABLE OF CONTENTS

Revision History	3
Damage in Transportation	4
Language Policy	5
Section 1 - Introduction	13
1-1 Product Identification and Shipping List.....	13
1-2 Compatibility	13
1-2-1 System Compatibility Warning.....	14
1-3 Related Documentation	14
1-4 Environmental Requirements	14
1-5 Theory of Operation.....	14
Section 2 - Setup and Calibration	16
2-1 Coil Installation	16
2-1-1 Special Install Notes	16
2-1-2 Installing the Coil	16
2-2 Installation Functional Checks.....	16
2-3 Periodic Quality Assurance Check	16
Section 3 - Functional Checks	17
3-1 Scanner Verification	17
3-2 Coil Imaging Performance Verification	17
3-2-1 Tools Required	17
3-2-2 Explanation of Procedure	17
3-2-3 Signal Scan.....	17
3-2-4 Noise Scan	20
3-2-5 SNR Image Analysis.....	21
3-3 External Cable Check.....	22
3-4 PIN Diodes Check	23
Section 4 - Maintenance	24
4-1 Coil Care.....	24
4-2 Special Care Requirements (if applicable)	24
Section 5 - Replacement.....	25
5-1 Assembly/ Disassembly of Coil.....	25
5-2 External Cable Replacement.....	26
5-3 Mechanical Hardware Replacement.....	26
Section 6 - Renewal Parts	28
6-1 Field Replaceable Units.....	28
Section 7 - Appendix.....	29
7-1 SNR Data Sheet	29
7-2 Coil Configuration	30

REVISION HISTORY

REV	DATE	AUTHOR	REASON FOR CHANGE
0	Mar. 28, 2002	G. Shajan	Initial Release
1	Apr. 12, 2002	G. Shajan	GE Catalogue Part No. changed from M1089GP to M1085GF (Sheet 1)
2	May 7, 2002	G. Shajan	Changed the header reference; Updated system references on front page; Updated system compatibility in Section 1-2; Updated CRM Configuration in Section 7-2.
3	Nov. 5, 2002	G. Shajan	Added Coil Non-Compatibility with Signa Advantage Systems with warning symbol on Page 7.
4	Oct. 15, 2007	Lalit Rai	Removed Service Manual from shipping list Table 1-1; Added HDx and HDe 1.5T Compatibility in Section 1.2; Updated O/M P/N in Section 1.3; Updated weight of Coil in Section 1.4; Updated Dimensional Specification with Strap Length and corrected Length from 56 cm to 55.6 cm
5	Nov. 15, 2007	Lalit Rai	Dimensions corrected in Section 1-1.
6	May 28, 2010	R. David	Removed PIN Diode Replacement section and P/N Diode FRU per SPRs MRIhc47186 and MRIhc44718.

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, the notation "Damage in Shipment" should be written on all copies of the freight or express bill before delivery is accepted or signed for by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14-day period.

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ENGLISH

(EN)

THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.

- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

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
SECTION 1 - INTRODUCTION

1-1 Product Identification and Shipping List

The coil consists of one major assembly (Illustration 1-1). The unit is 55.6 cm long (not including straps) by 22 cm wide. It has two mesh "windows" that indicate the open center of the two coils and aid in positioning.

Two 78.5 cm long straps are attached to one end of the coil fabric. The underside of these straps is covered with Velcro®, which attach to a Velcro hook strip at the other end of the coil and to four Velcro hooks located above and below the mesh windows.

A 183.3 cm long RF cable exits the material on one of the long sides and has a BNC connection that attaches to the Surface Coil quick-disconnect box.

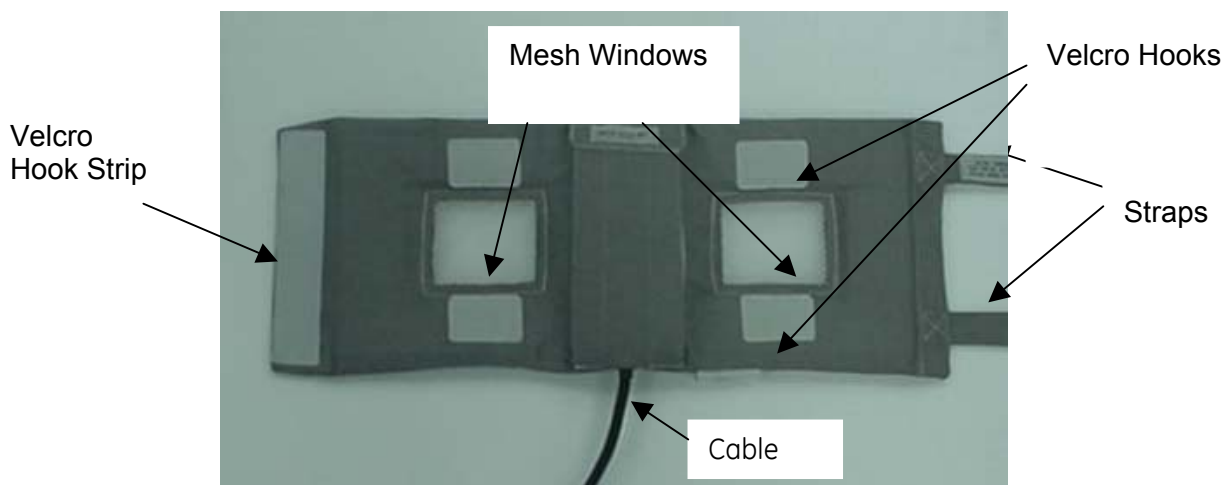GP FLEX COIL
ILLUSTRATION 1-1

TABLE 1-1
SHIPPING LIST

Description	GE Part#	Qty
1.5T GP Flex Coil	2320288	1
Operator Manual-CD	2336151-247	1

1-2 Compatibility

This coil is compatible with the following hardware configurations:

- Signa Horizon 1.5T
- Signa Horizon LX 1.5T
- Signa Twin Speed 1.5T
- Signa Infinity 1.5T
- Signa Infinity Twin Speed 1.5T
- Signa CV/i 1.5T
- Signa MR/i 1.5T
- Signa HDx 1.5T
- Signa HDe 1.5T

1-2-1 System Compatibility Warning

Do not use this coil on Signal Advantage 1.5T Systems. Patient burns may result.

1-3 Related Documentation

- Signal LX Service Methods CD, 2160623-1
- 1.5T GP Flex Coil Operator Manual, 2336151-247

1-4 Environmental Requirements

Storage Requirements

Coil should be stored in the Scanner Room.

Dimensions

Dimensions (excluding cable assembly): 55.6 cm x 22 cm x 4 cm, Strap Length: 78.5 cm.

Weight

Coil Weight: 575 grams maximum


1-5 Theory of Operation

The Signal General Purpose (GP) Flex coil design is a linear, receive-only flexible coil. The GP Flex Coil consists of two loops that are serially connected to create a co-rotating "saddle coil" pair. When aligned properly, this "saddle" design provides exceptional uniformity across the sensitive volume of the coil and a minimum of bright spots or high-intensity areas. If the loops are positioned in a more open fashion (i.e., they move farther away from each other, toward a planar, or flat configuration), the image uniformity decreases. For close positioning do not overlap the coil ends. Doing so will not damage the coil, but the performance will suffer. For optimal positioning the circular shape should be maintained with the two-mesh windows 10 cm apart and the coil ends 5 cm apart. Some deviation from the optimal position will not significantly degrade the uniformity. Therefore, this coil can be used in anatomical areas such as the hip, brachial plexus, larger knees and ankle.

The coil consists of one major assembly (see Illustration 1-1). The two loops are in a flexible material, covered with fabric. This allows the GP Flex Coil to bend in one direction, making it easy to wrap around the anatomy of interest.

The Functional Block Diagram of the antenna electronics is shown in Illustration 1-2. The two flex PCB wings on either side of the rigid PCB form the loops. The rigid PCB has a matching network that matches the loaded antenna to 50 ohms, and an active blocking network.

The blocking network consists of a single PIN diode and a tuned LC network. The diode is forward biased whenever the body coil transmits. The external cable is connected to the antenna via a right angle SMA jack. The flex PCB has a passive blocking network, which consists of a Schottky Diode, and a tuned LC network.

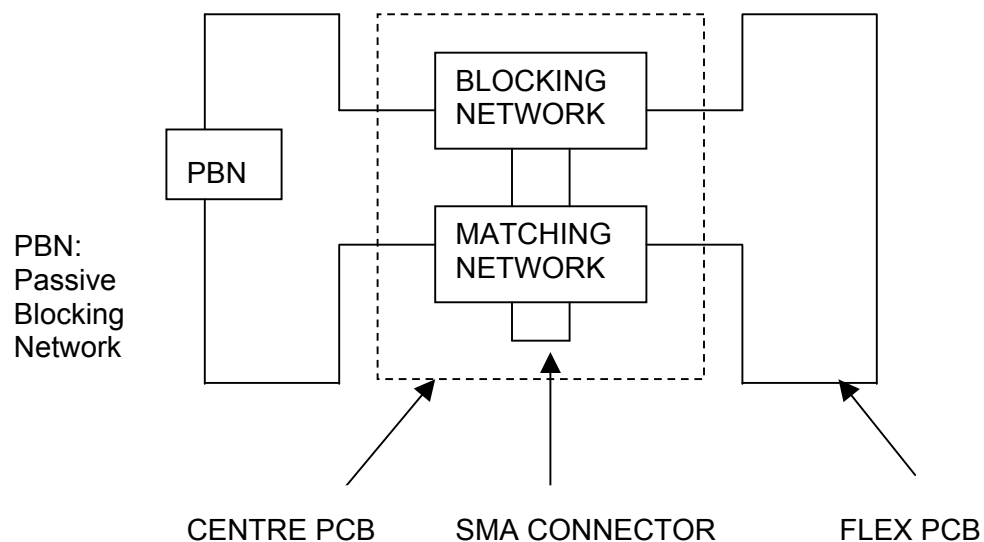GP FLEX COIL FUNCTIONAL BLOCK DIAGRAM
ILLUSTRATION 1-2

SECTION 2 - SETUP AND CALIBRATION

2-1 Coil Installation

2-1-1 Special Installation Notes

None

2-1-2 Installing the Coil

The name for this coil is: GPFLEX

Add the coil using the Configuration File Manager. Refer to the appropriate Service Methods CD > System Level Procedures > Software Utilities.

If the coil does not exist in Coil Config File, refer to the Adding New Coils to Config File Manager procedure and use the coil configuration information in *Section 7-3* of this manual.

2-2 Installation Functional Checks

1. Perform system level Signal to Noise Check. Refer to the appropriate Service Methods CD > System Level Procedures > Functional Checks > Signal-to-Noise Check.
2. Perform *Section 3 Coil Imaging Performance Verification*.

2-3 Periodic Quality Assurance Check

On a periodic basis, such as during planned maintenance, perform the quality assurance checks as outlined below to ensure the coils is operating properly.

1. Check external cable for cracks or cuts.
2. Perform *Section 3 Coil Imaging Performance Verification* and record data values on Data Sheet.

SECTION 3 - FUNCTIONAL CHECKS

3-1 Scanner Verification

Perform system level Signal to Noise Check. Refer to the appropriate Service Methods CD > System Level Procedures > Functional Checks > Signal-to-Noise Check.

3-2 Coil Imaging Performance Verification

3-2-1 Tools Required

TABLE 3-2-1
TOOLS REQUIRED

Description	GE Part#	Qty
Head TLT Sphere	46-265826G6	1
Head Loader	46-287899G1	1

3-2-2 Explanation of Procedure

SNR measurements for the GPFLEX coil require sets of signal and noise scans. Refer to the *Data Sheet* in Appendix 7-1 to understand the data required to calculate the SNR. All ROI measurements are made on the individual element images, **not** on the composite image.

The image quality check uses two different protocols for signal and noise image acquisition. The signal scan is an **FSE** sequence used to minimize susceptibility and B₀ inhomogeneity effects. The noise scan is a **GRE** sequence that has a Control Variable (do_noise) to eliminate the transmit RF completely during the scan. The signal scan **must** be run prior to the noise scan as the R1, R2, and TG values from the signal scan are used for the noise scan.

3-2-3 Signal Scan

The following procedure is specific to the LX platform but can be easily adapted for 5.x systems.

1. The Quad Head Coil must be completely removed from the cradle before performing any body or surface coil scans. Failure to do this may result in damage to the Quad Head Coil T/R network.
2. At the magnet, press the **[Alignment Light]** button to turn on the light. Move the cradle to align the coil to the alignment lights as shown in illustration 3-1. Press the **[Landmark]** button to landmark the alignment.
3. Move the coil to scan position by pushing the **[Move to Scan]** button, ensuring cable does not get snagged.
4. At the console, set the protocols as listed below.

SIGNAL PROTOCOL

Patient/Exam. Information

ID: **geservice**
Name: **GP Flex coil**
Patient Weight: **111 lbs. (50 kg)**

Patient Position

Patient position: **Supine**
Patient entry: **Head First**
Coil: **GP Flex**
Series Description: **Signal**

Imaging Parameters

Plane: **Axial**
Mode: **2D**
Pulse Seq.: **FSE**
Imaging Options: **Fast**
PSD Name: *leave blank*
Protocol: *leave blank*

Scan Timing

of Echoes: **1**
TE: **17**
TR: **500**
Echo Train Length: **4**
Bandwidth: **15.63**

Additional Parameters

No entries required in this area.

Acquistion Timing

Freq.: **256**
Phase: **256**
NEX: **1**
Phase FOV: **1.00**
Freq. DIR: **R/L**
Auto Center Freq.: **Peak**
Autoshim: **On**
Phase Correct: **On**
Contrast: **Off**
of Reps B4 Pause: **0**

Scanning Range

FOV: **24**

Slice Thickness: **5**

Spacing: **32.5**

A/P Center R/L Center: **0.00**


S/I Start: **S75**

S/I End: **I75**

Slices: **5**

Table Delta: **0.00**

5. Click **[Save Series]** to download the protocols, then click **[Prepare to Scan]**.
6. Run **[Auto Prescan]**. Record the R1, R2 and TG values on the SNR Data Sheet (found at the end of this manual).
7. Run **[Scan]**.

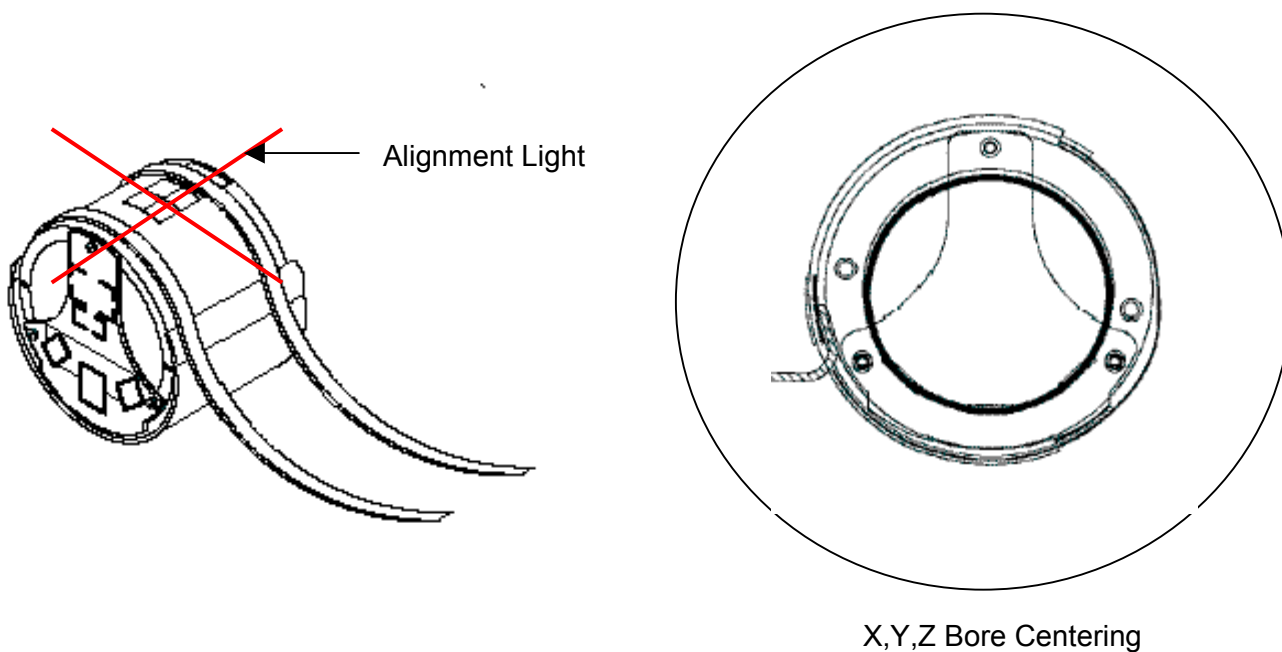LANDMARK AND ALIGNMENT
ILLUSTRATION 3-1

3-2-4 Noise Scan

A signal scan must be run **prior** to the noise scan as the same R1, R2 and TG values must be used for both the signal and noise scans. Do **not** run an Auto Prescan prior to the Noise Scan as the values will be changed.

1. Copy the signal scan series. Use **[Copy Series]** (highlight signal series and right mouse click the button) and **[Paste Series]** in RX Manager.
2. Click **[View Edit]** and set the protocols as listed below.

Patient Position

Patient position: **Supine**
Patient entry: **Head First**
Coil: **GP Flex**
Series Description: **Noise**

Imaging Parameters

Plane: **Axial**
Mode: **2D**
Pulse Seq: **GRE**
Imaging Options: **None**
PSD Name: *leave blank*
Protocol: *leave blank*

Scan Timing

of Echoes: **1**
TE: **Min Full**
TR: **34**
Flip Angle: **1**
Bandwidth: **15.63**

Additional Parameters

No entries required in this area.

Acquisition Timing

Freq: **256**
Phase: **256**
NEX: **1**

Phase FOV: **1.00**
Freq DIR: **R/L**
Auto Center Freq: **Peak**
Autoshim: **On**

Scanning Range

FOV: **24**

Slice Thickness: **5**

Spacing: **32.5**

A/P Center R/L Center: **0.00**

S/I Start: **S75**

S/I End: **I75**

Slices **5**

Table Delta: **0.00**

3. Click **[Save Series]** and click **[Prepare to Scan]**.
4. Open *Display CVs* menu under Research Operations. Set the rhformat and do_noise CVs to **1**.
5. Run **[Manual Prescan]**, do **not** make any changes, and click **[Done]**.
6. Run **[Scan]**.

3-2-5 SNR Image Analysis


SNR Measurement

Regions of interest in both signal and noise images can be measured directly in the image browser. Click the user interface **[Measure]** button, select the circular or rectangular shape, and adjust its size and orientation when the shape is displayed in the selected image. Mean, standard deviation, and area of the ROI will display in the lower right corner of the image.

$$\text{SNR}_{\text{Element}} = \frac{\text{Mean of Signal Within ROI}}{\text{Std. Dev. of Noise Within ROI}}$$


Note: The SNR calculation uses the MEAN of the signal image and STANDARD DEVIATION of the noise image. SNR is measured for **each element**, not on the composite image.

For the signal measurement, choose an ROI covering approximately 80% of the phantom. Also, measure the noise with an ROI covering 80% of the FOV. Examples of typical ROIs are shown in illustration 3-2 and 3-3.

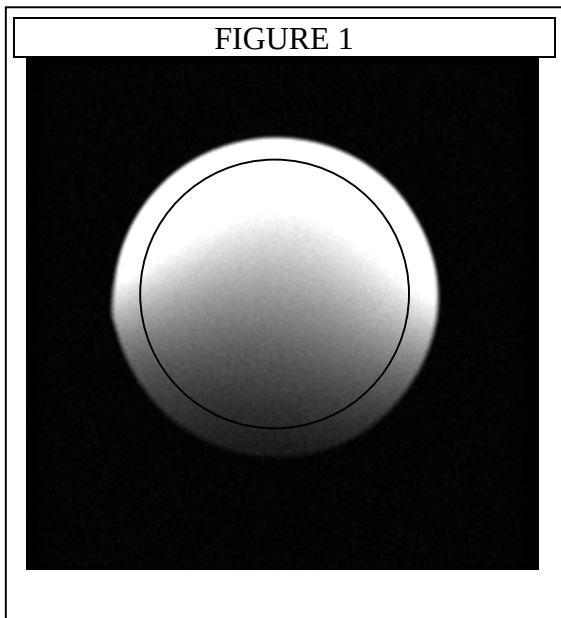SIGNAL SCAN

ILLUSTRATION 3-2

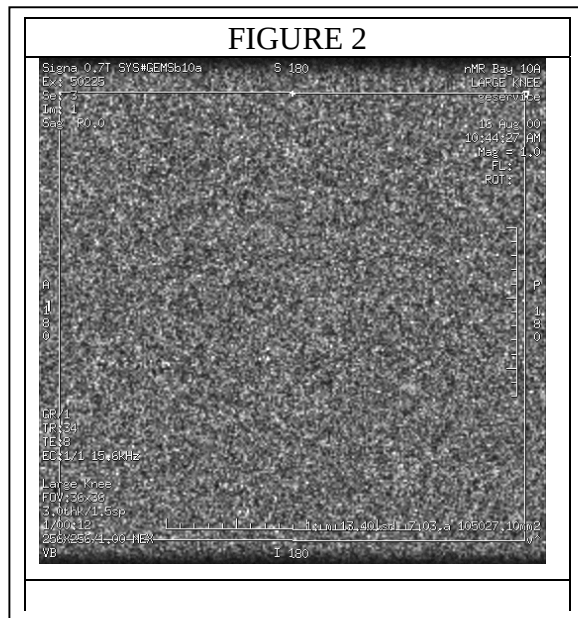NOISE SCAN

ILLUSTRATION 3-3

SNR Specification

The SNR measurements must be greater than or equal to the following specifications:

TABLE 3-2-5
SNR SPECIFICATIONS


Slice Location	Image Number	SNR
S37.5	2	72.5
S0	3	106
I37.5	4	72

Note: Image 1 and image 5 are collected to verify that the phantom is correctly centered. For example, images 1 and 5 should be the same size. If they are not the same size, then the phantom should be repositioned relative to the coil center.

3-3 External Cable Check

1. Select the DIODE TEST function on the DMM.
2. Connect the POSITIVE lead of the DMM to the center pin of the BNC on the external cable. Connect the NEGATIVE lead to the main housing/ground of the BNC.
3. Flex the external cable, especially near the connectors and the strain relief, and observe that a reading of 0.400 to 0.600 should remain on the DMM with no instability or fluctuations.

4. Connect the POSITIVE lead of the DMM to the main housing/ground of the BNC on the external cable. Connect the NEGATIVE lead to the center pin of the BNC.
5. Flex the external cable, especially near the connectors and the strain relief, and observe that a reading of INFINITY should remain on the DMM with no instability or fluctuations.

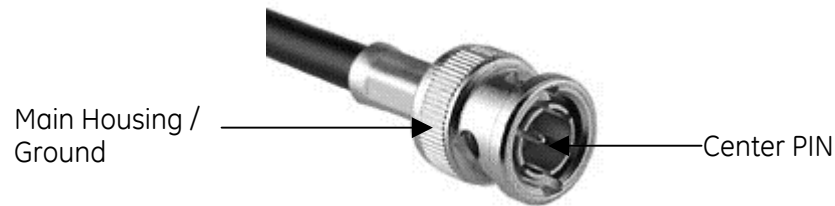BNC CONNECTOR
ILLUSTRATION 3-4

6. If the cable fails any of the above tests, replace it. Refer to *Section 5-2 External Cable Replacement*.

3-4 PIN Diode Check

Note: There is only one PIN diode in this antenna. This procedure will indicate if the PIN diode is defective.

1. Select the DIODE TEST function on the Digital Multi-meter (DMM).
2. Connect the POSITIVE lead of the DMM to the center pin of the BNC on the external cable. Connect the NEGATIVE lead to the main housing/ground of the BNC.
3. A reading of 0.400 to 0.600 should be observed on the DMM.
4. If a reading below 0.400 is observed in either direction, either the output cable is shorted or the PIN diode is bad.
5. If a reading above 0.600 is observed in Step 3, the PIN diode is defective.
6. Connect the POSITIVE lead of the DMM to the main housing/ground of the BNC on the external cable. Connect the NEGATIVE lead to the center pin of the BNC.
7. A reading of INFINITY should be observed on the DMM.
8. If a reading of INFINITY is observed in both directions, either the output cable is open or the PIN diode is open.
9. If any of the above conditions fails, the PIN diode is faulty. Due to safety hazards arising from soldering in the field, PIN diodes are no longer FRUs. If the coil has a defective PIN diode, the respective coil component must be replaced.

SECTION 4 - MAINTENANCE

4-1 Coil Care

WARNING!

Detach coil connector from scanner before attempting to clean. Do not reattach after cleaning until coil has dried completely. Having the coil attached to the system during cleaning or when it is wet may result in electrical shock.

CAUTION

Do not spray or pour cleaning solution directly on coil. Do not submerge coil in solution. The coil contains sensitive electronics components that could be damaged by the solution.

If necessary, a 10% bleach solution can be used for more rigorous cleaning; however, some discoloration of the strap, cover or pad may occur.

4-2 Special Care Requirements (If Applicable)


Contact the GE authorized service personnel for disposing the coil at the end of product life.

SECTION 5 - REPLACEMENT

Simple removals that are clearly obvious are not described here.

Unless otherwise noted, the steps for reassembly are simply the reverse order of the steps described for disassembly.

5-1 Assembly/ Disassembly of Coil

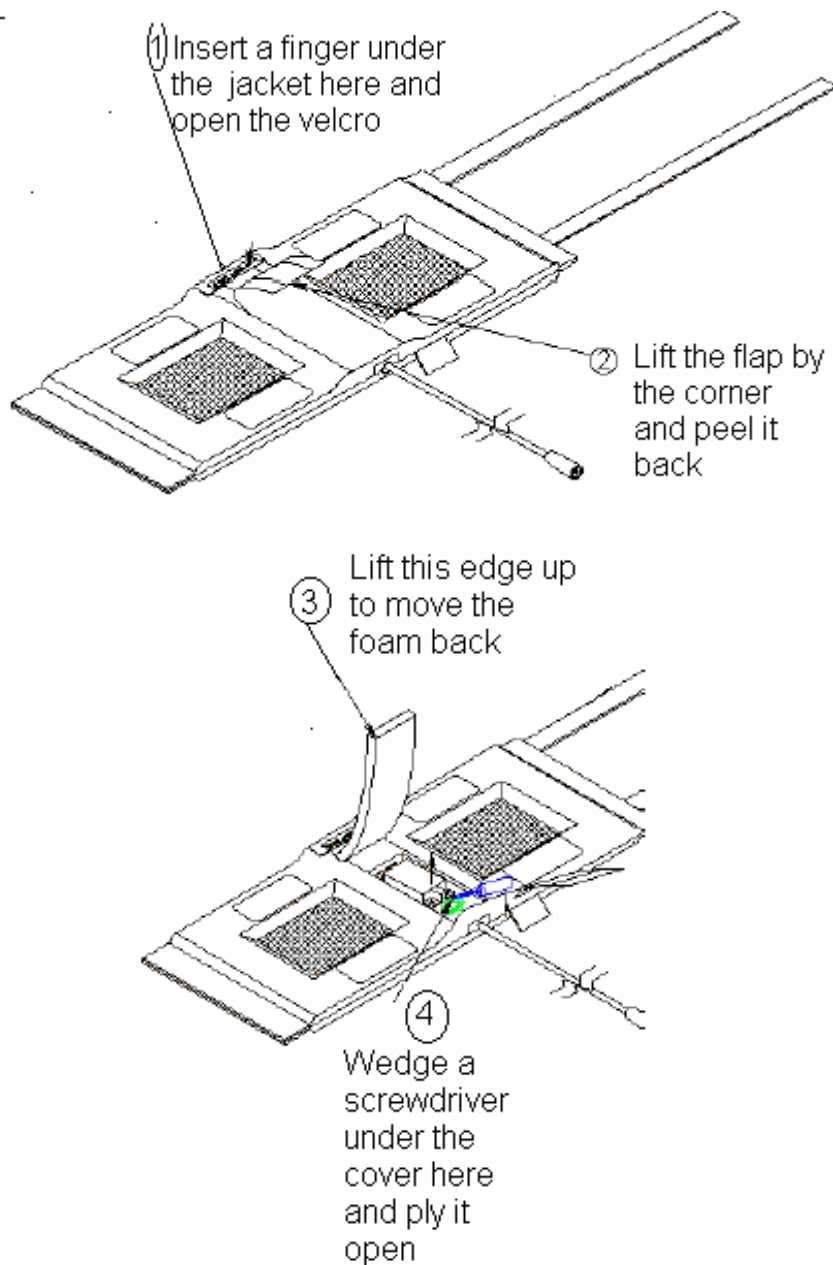
COIL ASSEMBLY / DISASSEMBLY
ILLUSTRATION 5-1

To reassemble, replace the cover on the posts, and press on it until it snaps into place. Push the foam flap under the jacket. Carefully put the fabric flap into place, making sure that the Velcro on the edges lines up with the Velcro on the jacket. Tuck the end of the flap under the jacket. Press over the Velcro while making sure that **no** edges are left open.

5-2 External Cable Replacement

1. Disassemble the coil as per *Section 5-1 Disassembly/ Assembly of Coil*. Unscrew the two nylon screws and remove the strain relief upper cover.
2. Unscrew (use two 8 mm wrenches, one on the circuit board female connector and one on the male cable connector) and remove the male SMA cable connector from the circuit board. Be careful not to twist the female SMA connector from the circuit board.
3. Put the external cable out of the coil.
4. Install the new cable. Stabilize the circuit board and the connector on the circuit board. Too much torque (>2 kgf-cm) can cause damage.
5. Put the strain relief back on. Reassemble the coils as per *Section 4-1 Disassembly/ Assembly of GP Flex Coil*.

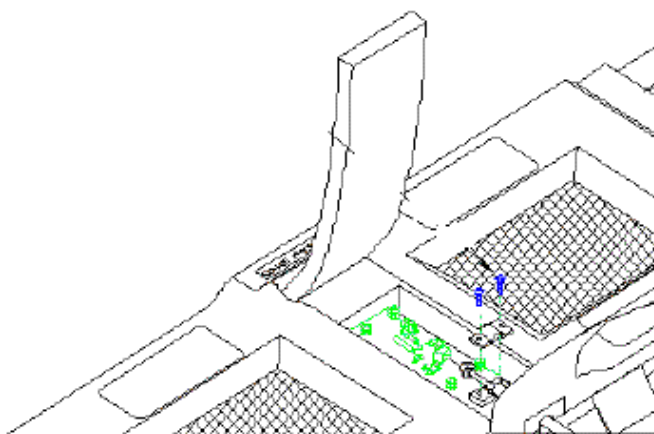
REPLACING EXTERNAL CABLE
ILLUSTRATION 5-2

5-3 Mechanical Hardware Replacement

1. Refer to *Section 6 Renewal Parts* for the hardware kit part number.
2. To replace the cover assembly refer to *Section 5-1 Disassembly/Assembly of Coil*.
3. To replace the upper strain relief block and the two screws, refer to *Section 5-2 Replacing the External Cable*.

4. If a plastic standoff breaks, it can be pulled out using a needle nose pliers. Insert a new standoff in its place, making sure that the flat end goes in first. See Illustration 5-3.

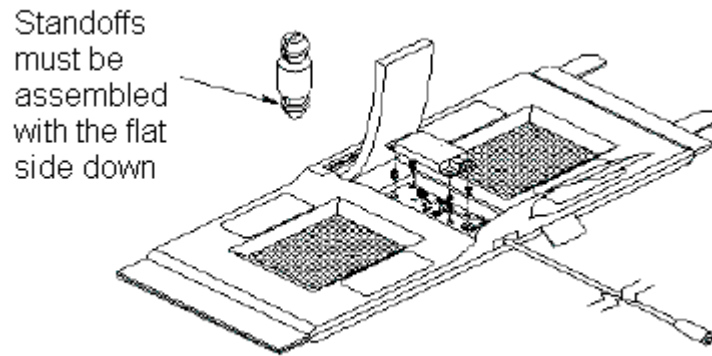REPLACING MECHANICAL HARDWARE
ILLUSTRATION 5-3


SECTION 6 - RENEWAL PARTS

6-1 Field Replaceable Units


TABLE 6-1
FIELD REPLACEABLE UNITS LIST

Description	GE Part#
1.5T GP Flex Coil	2320288
CABLE KIT	2337313
HARDWARE KIT	2337313-2


1.5T SIGNA GP FLEX COIL WITH CABLE, P/N 2320288

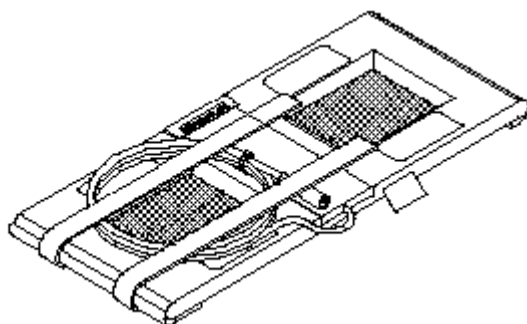1.5T SIGNA GP FLEX COIL CABLE KIT, P/N 2337313

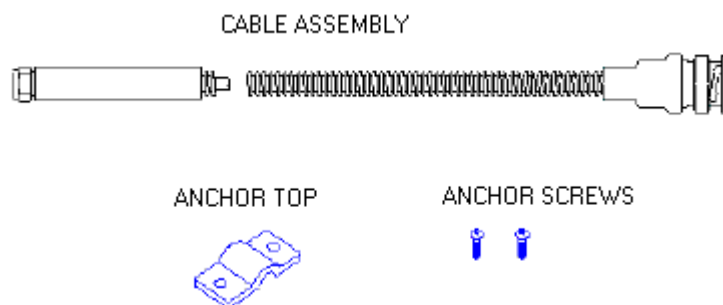1.5T SIGNA GP FLEX COIL HARDWARE KIT, P/N 2337313-2

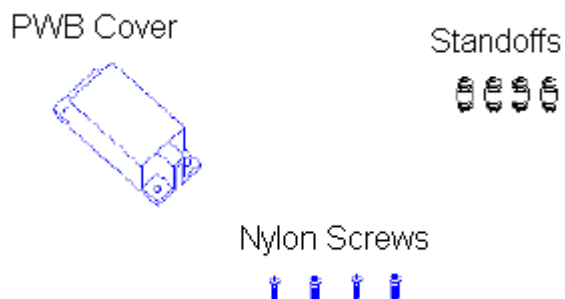SECTION 7 - APPENDIX

7-1 SNR Data Sheet

Use the table provided below to record the calculated signal to noise ratio (SNR) data obtained from the Functional Checks section.

GP Flex Coil SNR Datasheet

Date	Comments						
Slice Location	R1	R2	TG	Signal Mean (A)	Noise Std. Dev. (B)	SNR (A/B)	SNR Spec Limit
S37.5							72.5
S0							106
I37.5							72

Date	Comments						
Slice Location	R1	R2	TG	Signal Mean (A)	Noise Std. Dev. (B)	SNR (A/B)	SNR Spec Limit
S37.5							72.5
S0							106
I37.5							72

Date	Comments						
Slice Location	R1	R2	TG	Signal Mean (A)	Noise Std. Dev. (B)	SNR (A/B)	SNR Spec Limit
S37.5							72.5
S0							106
I37.5							72

7-2 Coil Configuration

Parameter	1.5T CRM Configuration	1.5T BRM Configuration
Coil Name	GPFLEX	GPFLEX
Coil Type	3	3
Extremity Coil	No	No
Cable Loss	1.2	1.05
Coil Loss	0.952	1.72
Recon Scale Factor	0.75	0.75
Linear vs. Quadrature	1	1
Multiple Receiver Coil	No	No
Number of Receivers	0	0
Starting Receiver ID	0	0
Ending Receiver ID	0	0
Multi-Coil Port Enable	0	0
Multi-Coil Part Error Enable	0	0
Additional Transmit Attenuation	0	0
Number of Fast Receivers	1	1
Starting Fast Receiver ID	4	4
Ending Fast Receiver ID	4	4
Start TA Value	90	90
Start RG Value	12	12
Multi-Coil Recon Enable	0	0
Head Default Freq. Direction	0	0
SCIC Axial	0	0
SCIC Sagittal	0	0
SCIC Coronal	0	0
Auto Shim Receiver	-1	-1

Page intentionally left blank.

www.gehealthcare.com

Imagination at work