

SECTION 6 – INTERCONNECT DATA

TABLE OF CONTENTS

<u>SECTION</u>	<u>TITLE</u>	<u>PAGE</u>
6-1	INTRODUCTION	6-3
6-1-1	Component Designators	6-3
6-1-2	Definition of Terms	6-4
6-2	SYSTEM INTERCONNECTS	6-5
6-2-1	Introduction	6-5
6-2-2	M1085KF – Mobile Water-to-Water Heat Exchanger Kit	6-6
6-2-3	M1090YL – Signa MR/i LCC Magnet WideOpen Enclosure Upgrade Site Collector Kit	6-6
6-2-4	M1090YN – Signa MR/i True Mobile LCC Magnet WideOpen Enclosure Upgrade Site Collector Kit	6-6
6-2-5	M1090YW – LCC Magnet WideOpen Enclosure Upgrade Kit	6-6
6-2-6	M1091YM – Signa CV/i LCC Magnet WideOpen Enclosure Upgrade Site Collector Kit	6-7

6-1 INTRODUCTION

6-1-1 Component Designators

GE uses a Component Designator System as a means of identifying system components in a consistent manner. All subsystem cabinets and other components are referred to by their component designators in the diagrams and tables of this section. For example, the System Cabinet is referred to as MR2. Refer to Table 6-1 for all component designators.

TABLE 6-1
COMPONENT DESIGNATION

BASIC SYSTEM OR OPTION	COMPONENT DESIGNATOR	DESCRIPTION
System Options	EO1/EO2	Emergency Off Buttons
	HE1	Mobile System Water-to-Water Heat Exchanger
	MG2	Magnet Enclosure
	MG3	Magnet Rear Pedestal
	MG6	Blower Box
	MR1	RF/Penetration Cabinet for 1.5T system RF/PDU Cabinet for 1.0T system (PD1 is a module in lower portion of cabinet)
	MR2	System Cabinet
	MR3	SGD Base Gradient Cabinet for Base System or SGD Hi-Slew Gradient Cabinet for HiSpeed or EchoSpeed System
	MS1	Superconducting Magnet
	MS4	Magnet Rundown Unit
	MS5	Shield/Cryo Cooler Compressor Cabinet
	MSM	Magnet Monitor
	OW1	Operator Workspace
	PD1	Compact Power Distribution Unit for 1.5T system ONLY
	PA1	Pneumatic Patient Alert Control Box
	PP1	Penetration Panel
	PT1	Patient Transport Table
	WC1	Water Chiller for Gradient Coil (For Fixed Sites Only)
System Options	OM1	Oxygen Monitor
	OM3	Remote Oxygen Sensor Module
System CV/i Options	ADW1	LCD Color Monitor
	ADW2	Keyboard
	ADW3	Ultra2 Sparc Computer

6-1-2 Definition of Terms

The definition of terms used in following sections tables are:

- Group Number: identifying number referenced to bundles (i.e. groups) of cables as shown in interconnect diagrams
- Between Units (From/To): component designators as found in tables
- Usable Length: total length of a cable MINUS any required take up within cabinets
- Run Number: unique number assigned to each GE-supplied cable


Note

The run number must be used when making special cable orders.

- Cable Diameter: diameter of an individual cable
- Plug Pulling Diameter x Length: cable plug dimensions

Note

In some cases, a cable has more than one connector on an end. These cables will have the number of connectors following the english dimensions of the plug pulling diameter times length (e.g. '2.0x3.25 -x2' means there are 2 connectors with dimensions of 2.0 in. diameter and 3.25 in. length). Of course, the same number of connectors apply to the metric dimensions as well.

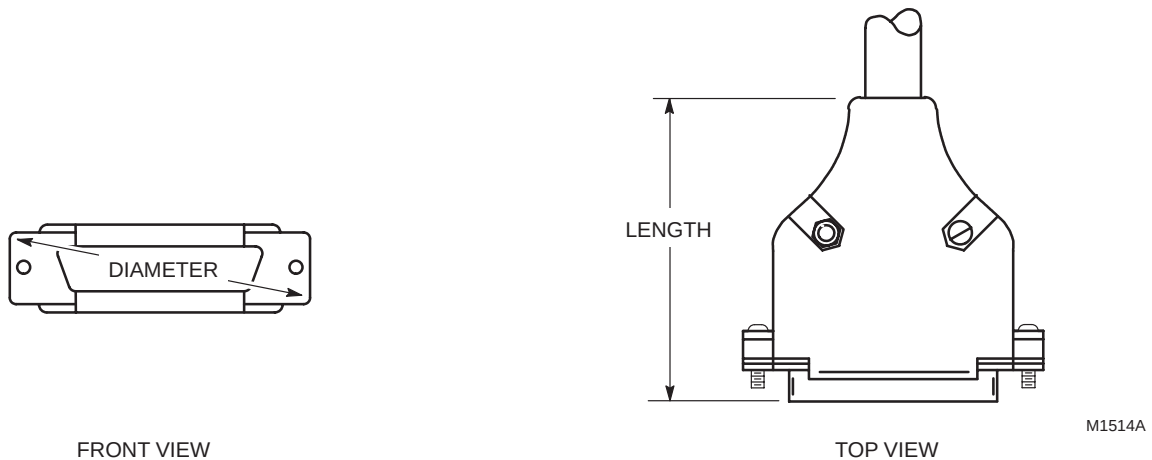


SUBMINIATURE-D CONNECTOR PLUG PULLING DIMENSIONS
ILLUSTRATION 6-1


6-2 SYSTEM INTERCONNECTS ADDED BY UPGRADE CATALOGS

6-2-1 Introduction

Illustrations 6-2 and 6-3 interconnect diagrams (1.5T systems and 1.0T systems) show the groups to which cables are added by the Signa MR/i & CV/i Upgrades for system with LCC Magnet. The following sub-sections list the cables provided by the various upgrade catalogs. Conduit or pipe is not recommended for cable runs since the system uses many prefabricated cables with large connectors. Unless otherwise specified, cables and components listed in tables in this section are supplied by GE.

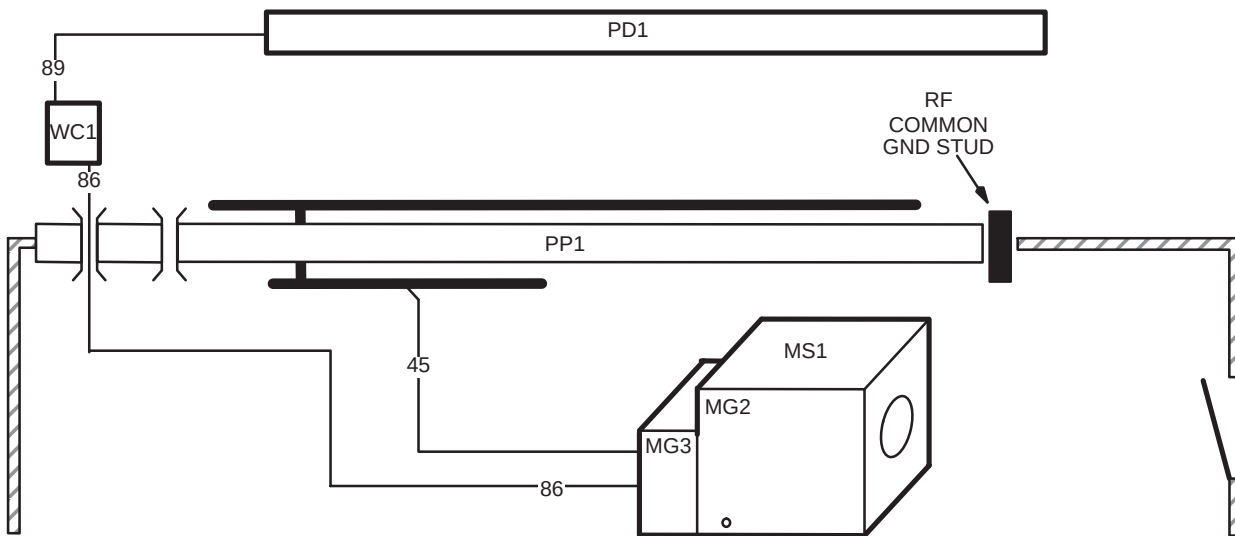

1.5T SYSTEM GROUP INTERCONNECT DIAGRAM
ILLUSTRATION 6-2

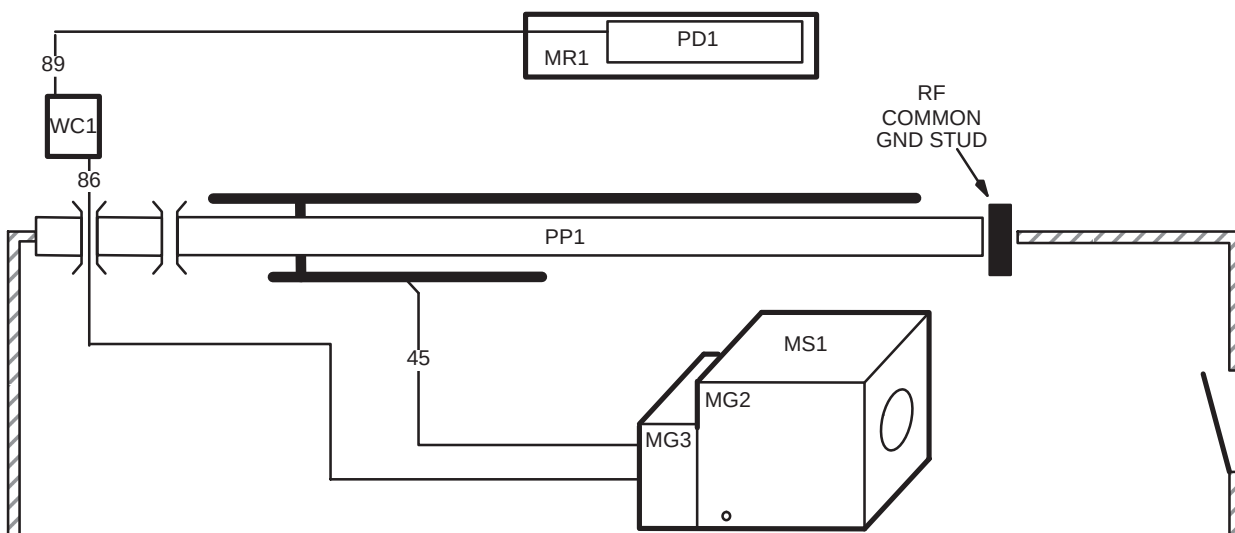

1.0T SYSTEM GROUP INTERCONNECT DIAGRAM
ILLUSTRATION 6-3

6-2-2 M1085KF – Mobile Water-to-Water Heat Exchanger Kit

The Mobile Water-to-Water Heat Exchanger Kit for Gradient Coil Cooling will be installed and connected to the van chilled water and the Signa system by the Van Manufacturer.

TABLE 6-2
CABLES PROVIDED BY CATALOG M1090YL

GROUP NUMBER	BETWEEN UNITS		USABLE LENGTH ft (m)	RUN NUMBER	CABLE DIAMETER in. (mm)	PLUG PULLING DIAMETER X LENGTH in. (mm)		NOTES
89	PD1	HE1	67 (20.3)	046	0.635 (16.1)	Hard Wired	Hard Wired	Mobile configurations do not use WC1.

6-2-3 M1090YL – Signa MR/i LCC Magnet WideOpen Enclosure Upgrade Site Collector Kit

TABLE 6-3
CABLES PROVIDED BY CATALOG M1090YL

GROUP NUMBER	BETWEEN UNITS		USABLE LENGTH ft (m)	RUN NUMBER	CABLE DIAMETER in. (mm)	PLUG PULLING DIAMETER X LENGTH in. (mm)		NOTES
45	PP1	MG3	30 (9.15)	300	0.425 (10.8)	1.30x2.00 (33.5x50.8)	Ring Terminals	
86	WC1	MG2	80 (24.4)	536 537	0.75 (19.1) 0.75 (19.1)	Flexible Tubing Flexible Tubing	Flexible Tubing Flexible Tubing	The existing system Runs 536 & 537 are re-used and must be rerouted to new Water Chiller for Gradient Cooling (WC1) from old (removed) Heat Exchanger (HE1).
89	PD1	WC1	67 (20.3)	065	0.635 (16.1)	Hard Wired	Hard Wired	Mobile configurations do not use WC1, van chiller system provides coil cooling water

6-2-4 M1090YN – Signa MR/i True Mobile LCC Magnet WideOpen Enclosure Upgrade Site Collector Kit

TABLE 6-4
CABLES PROVIDED BY CATALOG M1090YN

GROUP NUMBER	BETWEEN UNITS		USABLE LENGTH ft (m)	RUN NUMBER	CABLE DIAMETER in. (mm)	PLUG PULLING DIAMETER X LENGTH in. (mm)		NOTES
45	PP1	MG3	15 (4.56)	300	0.425 (10.8)	1.30x2.00 (33.5x50.8)	Ring Terminals	

6-2-5 M1090YW – LCC Magnet WideOpen Enclosure Upgrade Kit

TABLE 6-5
CABLES PROVIDED BY CATALOG M1090YW

GROUP NUMBER	BETWEEN UNITS		USABLE LENGTH ft (m)	RUN NUMBER	CABLE DIAMETER in. (mm)	PLUG PULLING DIAMETER X LENGTH in. (mm)		NOTES
45	PP1	MG3	60 (18.2)	385	0.525 (13.3)	2.80x2.00 (70.4x50.8)	2.80x2.00 (70.4x50.8)	

6-2-6 M1091YM – Signa CV/i LCC Magnet WideOpen Enclosure Upgrade Site Collector Kit

TABLE 6-6
CABLES PROVIDED BY CATALOG M1091YM

GROUP NUMBER	BETWEEN UNITS		USABLE LENGTH ft (m)	RUN NUMBER	CABLE DIAMETER in. (mm)	PLUG PULLING DIAMETER X LENGTH in. (mm)		NOTES
45	PP1	MG3	30 (9.15)	300	0.425 (10.8)	1.30x2.00 (33.5x50.8)	Ring Terminals	
86	WC1	MG2	80 (24.4)	536 537	0.75 (19.1) 0.75 (19.1)	Flexible Tubing Flexible Tubing	Flexible Tubing Flexible Tubing	The existing system Runs 536 & 537 are re-used and must be rerouted to new Water Chiller for Gradient Cooling (WC1) from old (removed) Heat Exchanger (HE1).
89	PD1	WC1	67 (20.3)	065	0.635 (16.1)	Hard Wired	Hard Wired	Mobile configurations do not use WC1, van chiller system provides coil cooling water