



GE Medical Systems

Technical Publications

Direction 2256047

Revision 0

Signa® Carriage Cover Upgrade to Low Profile Carriage Cover (LPCC)

Copyright© 1999 by General Electric Company

Operating Documentation

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "**damage in shipment**" written on **all** copies of the freight or express bill **before** delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

Call Traffic and Transportation, Milwaukee, WI (414) 785-5052/8*323-5052 **immediately** after damage is found. At this time be ready to supply name of carrier, delivery date, consignee name, freight or express bill number, item damaged and extent of damage.

Complete instructions regarding claim procedure are found in Section "S" of the Policy & Procedure Bulletins.

3/12/92



GE Medical Systems

GE Medical Systems: Telex 3797371
P.O. Box 414, Milwaukee, Wisconsin 53201 U.S.A.
(Asia, Pacific, Latin America, North America)

GE Medical Systems — Europe: Telex 261794
Shortlands, Hammersmith, London W6 8BX U.K.

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

AVERTISSEMENT

- CE MANUEL DE MAINTENANCE N'EST DISPONIBLE QU'EN ANGLAIS.
- SI LE TECHNICIEN DU CLIENT A BESOIN DE CE MANUEL DANS UNE AUTRE LANGUE QUE L'ANGLAIS, C'EST AU CLIENT QU'IL INCOMBE DE LE FAIRE TRADUIRE.
- NE PAS TENTER D'INTERVENTION SUR LES ÉQUIPEMENTS TANT QUE LE MANUEL SERVICE N'A PAS ÉTÉ CONSULTÉ ET COMPRIS.
- LE NON-RESPECT DE CET AVERTISSEMENT PEUT ENTRAÎNER CHEZ LE TECHNICIEN, L'OPÉRATEUR OU LE PATIENT DES BLESSURES DUES À DES DANGERS ÉLECTRIQUES, MÉCANIQUES OU AUTRES.

WARNUNG

- DIESES KUNDENDIENST-HANDBUCH EXISTIERT NUR IN ENGLISCHER SPRACHE.
- FALLS EIN FREMDER KUNDENDIENST EINE ANDERE SPRACHE BENÖTIGT, IST ES AUFGABE DES KUNDEN FÜR EINE ENTSPRECHENDE ÜBERSETZUNG ZU SORGEN.
- VERSUCHEN SIE NICHT, DAS GERÄT ZU REPARIEREN, BEVOR DIESES KUNDENDIENST-HANDBUCH NICHT ZU RATE GEZOGEN UND VERSTANDEN WURDE.
- WIRD DIESE WARNUNG NICHT BEACHTET, SO KANN ES ZU VERLETZUNGEN DES KUNDENDIENSTTECHNIKERS, DES BEDIENERS ODER DES PATIENTEN DURCH ELEKTRISCHE SCHLÄGE, MECHANISCHE ODER SONSTIGE GEFAHREN KOMMEN.

AVISO

- ESTE MANUAL DE SERVICIO SÓLO EXISTE EN INGLÉS.
- SI ALGÚN PROVEEDOR DE SERVICIOS AJENO A GEMS SOLICITA UN IDIOMA QUE NO SEA EL INGLÉS, ES RESPONSABILIDAD DEL CLIENTE OFRECER UN SERVICIO DE TRADUCCIÓN.
- NO SE DEBERÁ DAR SERVICIO TÉCNICO AL EQUIPO, SIN HABER CONSULTADO Y COMPRENDIDO ESTE MANUAL DE SERVICIO.
- LA NO OBSERVANCIA DEL PRESENTE AVISO PUEDE DAR LUGAR A QUE EL PROVEEDOR DE SERVICIOS, EL OPERADOR O EL PACIENTE SUFRAN LESIONES PROVOCADAS POR CAUSAS ELÉCTRICAS, MECÁNICAS O DE OTRA NATURALEZA.

ATENÇÃO

- ESTE MANUAL DE ASSISTÊNCIA TÉCNICA SÓ SE ENCONTRA DISPONÍVEL EM INGLÊS.
- SE QUALQUER OUTRO SERVIÇO DE ASSISTÊNCIA TÉCNICA, QUE NÃO A GEMS, SOLICITAR ESTES MANUAIS NOUTRO IDIOMA, É DA RESPONSABILIDADE DO CLIENTE FORNECER OS SERVIÇOS DE TRADUÇÃO.
- NÃO TENHA TENTADO REPARAR O EQUIPAMENTO SEM TER CONSULTADO E COMPREENDIDO ESTE MANUAL DE ASSISTÊNCIA TÉCNICA.
- O NÃO CUMPRIMENTO DESTA AVISO PODE POR EM PERIGO A SEGURANÇA DO TÉCNICO, OPERADOR OU PACIENTE DEVIDO A CHOQUES ELÉTRICOS, MECÂNICOS OU OUTROS.

AVVERTENZA

- IL PRESENTE MANUALE DI MANUTENZIONE È DISPONIBILE SOLTANTO IN INGLESE.
- SE UN ADDETTO ALLA MANUTENZIONE ESTERNO ALLA GEMS RICHIEDE IL MANUALE IN UNA LINGUA DIVERSA, IL CLIENTE È TENUTO A PROVVEDERE DIRETTAMENTE ALLA TRADUZIONE.
- SI PROCEDA ALLA MANUTENZIONE DELL'APPARECCHIATURA SOLO DOPO AVER CONSULTATO IL PRESENTE MANUALE ED AVERNE COMPRESO IL CONTENUTO.
- NON TENERE CONTO DELLA PRESENTE AVVERTENZA POTREBBE FAR COMPIERE OPERAZIONI DA CUI DERIVINO LESIONI ALL'ADDETTO ALLA MANUTENZIONE, ALL'UTILIZZATORE ED AL PAZIENTE PER FOLGORAZIONE ELETTRICA, PER URTI MECCANICI OD ALTRI RISCHI.

警告

- ・ このサービスマニュアルには英語版しかありません。
- ・ GEMS以外でサービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- ・ このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないで下さい。
- ・ この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

注意：

- 本维修手册仅存有英文本。
- 非 GEMS 公司的维修员要求非英文本的维修手册时，客户需自行负责翻译。
- 未详细阅读和完全了解本手册之前，不得进行维修。
- 忽略本注意事项会对维修员，操作员或病人造成触电，机械伤害或其他伤害。

REVISION HISTORY

REV	DATE	PRIMARY REASON FOR CHANGE
Rev A	.. Sep 24, 1999	.. Preliminary Release
Rev 0	... Oct 6, 1999	... Initial Release

LIST OF EFFECTIVE PAGES

PAGE	REV	PAGE	REV	PAGE	REV	PAGE	REV	PAGE	REV
Title Page		0							
Damage in Trans.		–							
Direction 2128126	... 0*								
A		0							
1 to 15		0							

* This revision number/letter corresponds to the indicated document's revision control system.

TABLE OF CONTENTS

SECTION	TITLE	PAGE
SECTION 1 – INTRODUCTION		
1–1	LOW PROFILE CARRIAGE COVER (LPCC) OPTIONS INSTALLED BY THIS MANUAL	1
SECTION 2 – 1.5T AND 1.0T NON-MULTICOIL LPCC UPGRADE		
2–1	NON-MULTICOIL UPGRADE PROCESS	3
2–2	NON-MULTICOIL 1.5T AND 1.0T LPCC ASSEMBLY	5
SECTION 3 – 1.5T AND 1.0T MULTICOIL LPCC UPGRADE		
3–1	MULTICOIL UPGRADE PROCESS	9
3–2	MULTICOIL 1.5T AND 1.0T LPCC ASSEMBLY	13

SECTION 1 – INTRODUCTION

1-1 1.5T AND 1.0T MULTICOIL COMPATIBILITY OPTIONS INSTALLED BY THIS MANUAL

This manual describes installation of the following options:

Catalog M1090WK, LPCC Upgrade for 1.5T Horizon without Multicoil SECTION 2

Catalog M1890WK, LPCC Upgrade for 1.0T Horizon without Multicoil SECTION 2

Catalog M1090WL, LPCC Upgrade for 1.5T Horizon with Multicoil SECTION 3


Catalog M1890WL, LPCC Upgrade for 1.0T Horizon with Multicoil SECTION 3

This page left blank intentionally.

SECTION 2 – NON-MULTICOIL LPCC UPGRADE


2-1 NON-MULTICOIL UPGRADE PROCESS

2A – TRANSFER EXISTING PREAMPS AND CONNECT TO BNC INPUT CABLES

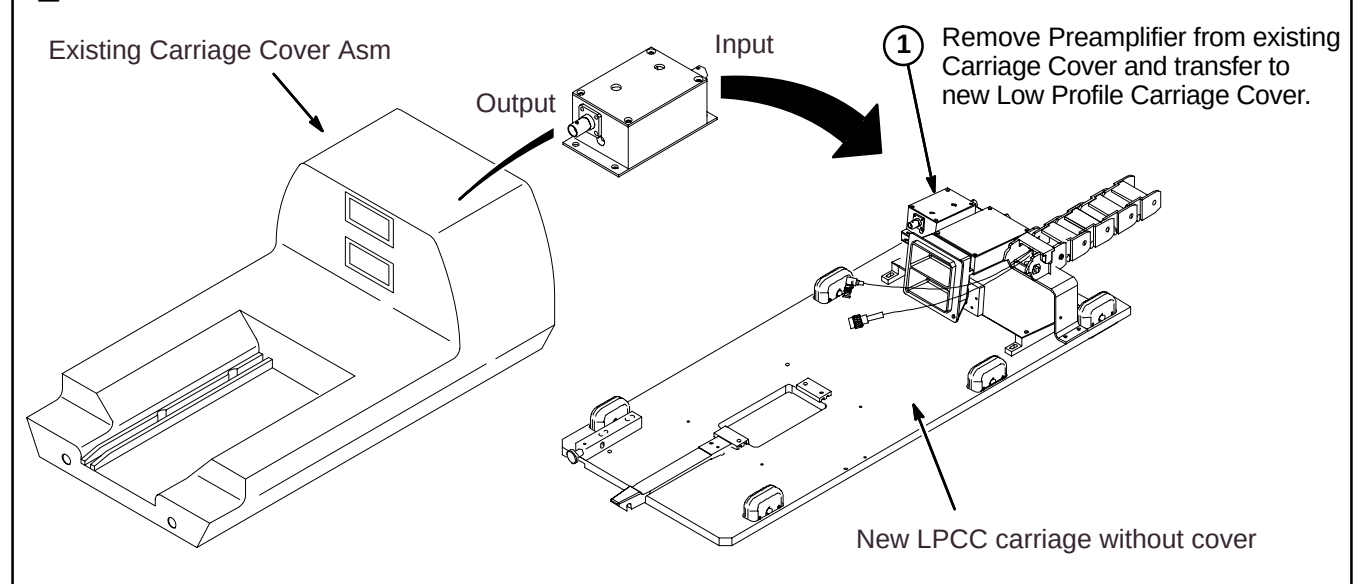


2B – CONNECT RUNS 404 AND 485

Note:
Cables must be pulled in or out of cable track.

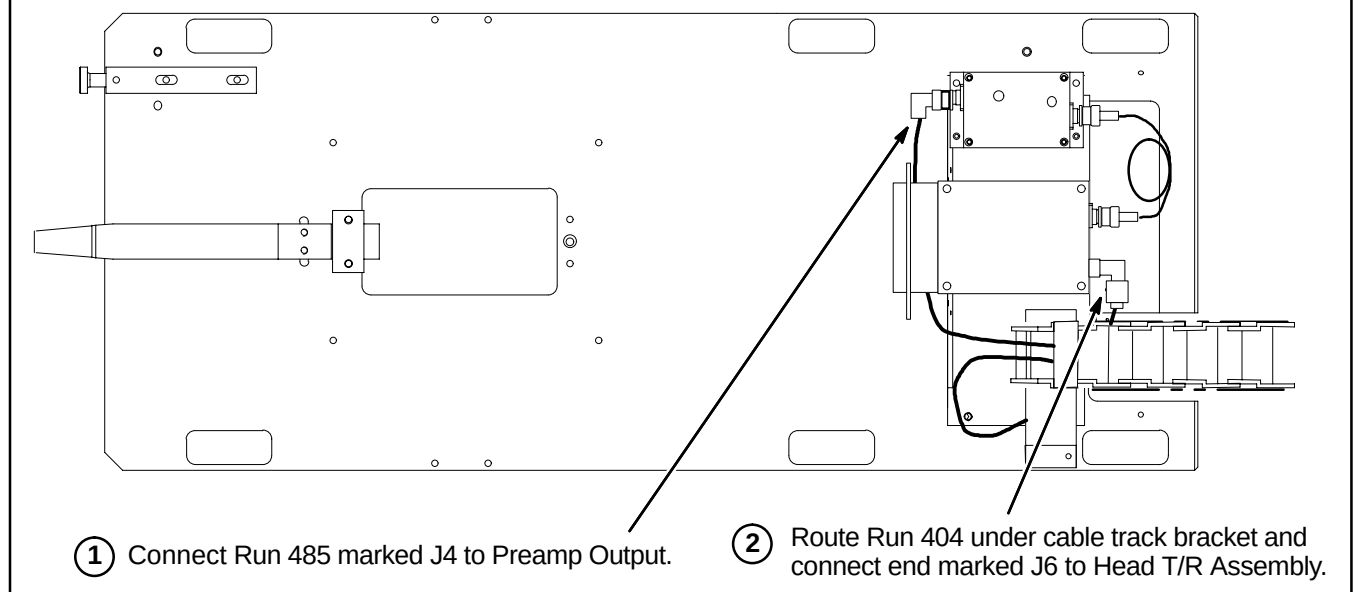

2-1 NON-MULTICOIL UPGRADE PROCESS (Continued)

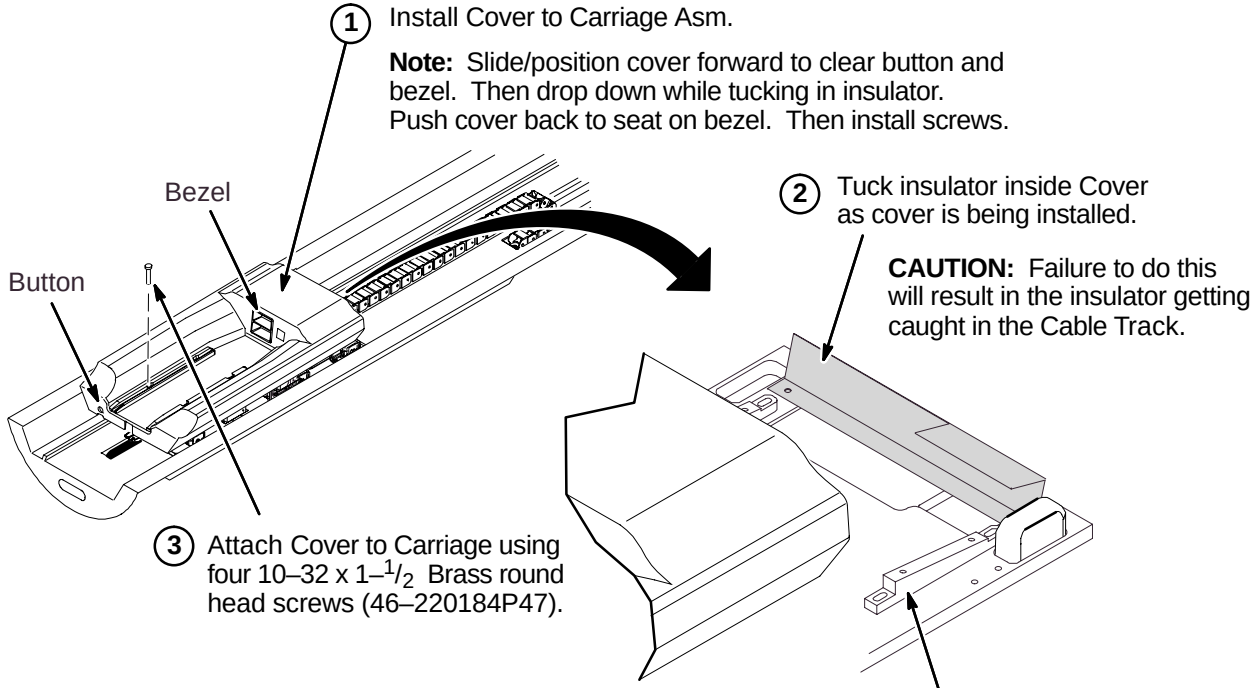A technical diagram showing the three-step process of installing a carriage cover. Step 1 shows the cover being positioned over the carriage assembly, with labels for 'Button' and 'Bezel'. Step 2 shows the cover being lowered, with an arrow indicating the direction of movement and a 'CAUTION' note about the insulator. Step 3 shows the cover being secured with screws. A large curved arrow connects the three steps. Labels include 'Button', 'Bezel', and 'Cable Track'.

2C – INSTALL CARRIAGE COVER

① Install Cover to Carriage Asm.

Note: Slide/position cover forward to clear button and bezel. Then drop down while tucking in insulator. Push cover back to seat on bezel. Then install screws.


② Tuck insulator inside Cover as cover is being installed.

CAUTION: Failure to do this will result in the insulator getting caught in the Cable Track.

③ Attach Cover to Carriage using four 10-32 x 1-1/2 Brass round head screws (46-220184P47).

Note: Mounting brackets may be adjusted forward or back for best fit of cover to base.

2-2 1.5T AND 1.0T NON-MULTICOIL LPCC ASSEMBLY

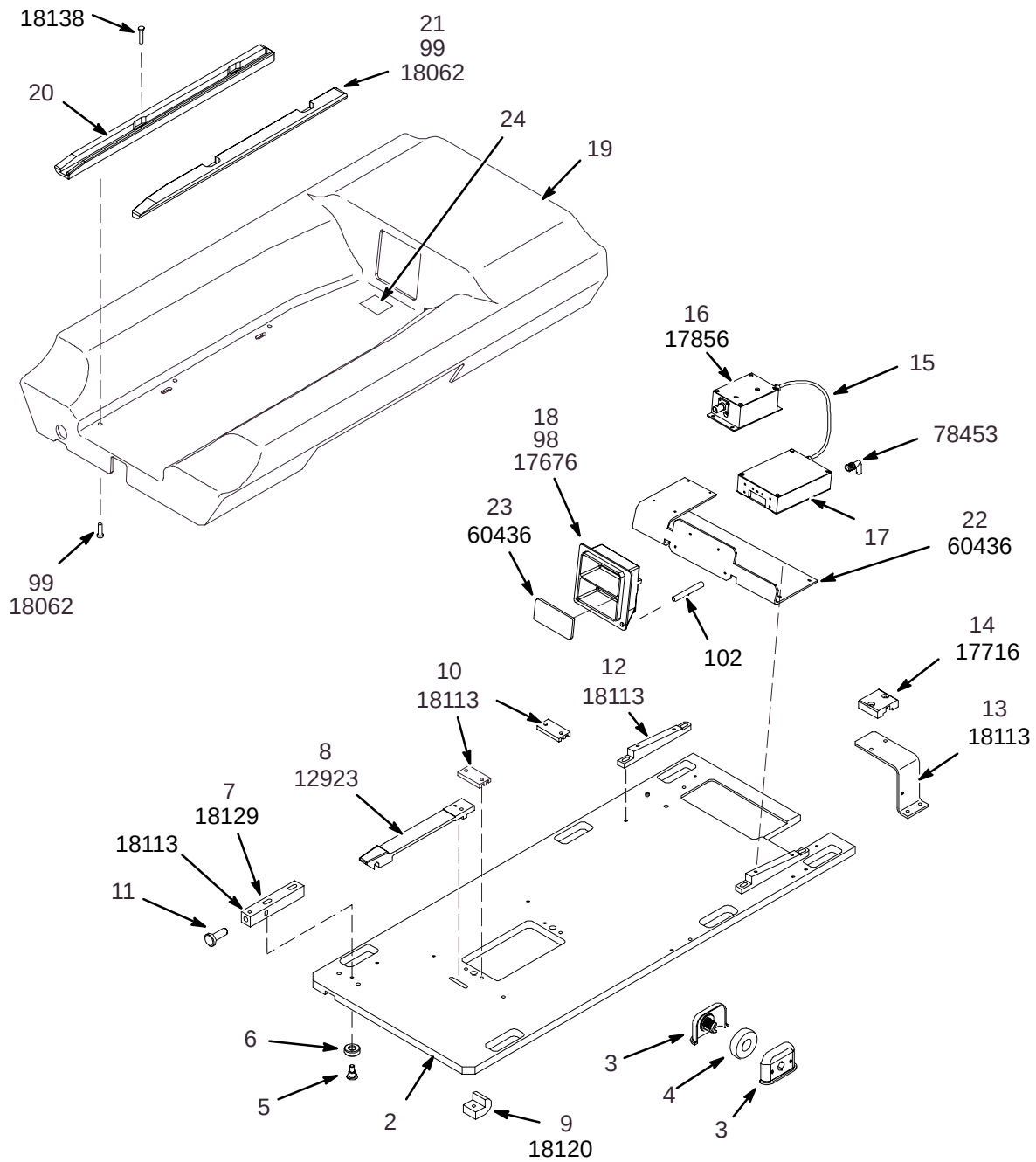

2-2 1.5T NON-MULTICOIL LPCC ASSEMBLY

Item	Part Number	FRU	Name	Quantity	Description (Remarks)
2	46-282163P2	2	BASE PLATE	1	LPCC CARRIAGE BASEPLATE
3	46-243167P1	2	BOX	12	WHEEL HOUSING, MOLDED NORYL
4	46-243494P1	2	CASTER	6	3/4" ID, 1-3/4 OD, 1/2" WIDE PLASTIC
5	46-243141P1	2	SHLD SCR	4	PLASTIC BEARING MOUNT SHOULDER SCREW
6	46-243178P2	2	BEARING	4	BEARING
7	46-282164P1	2	PLUNGER	1	CARRIAGE RELEASE PLUNGER, DELRIN
8	46-317169P1	2	HOOK	1	CRADLE LATCH HOOK
9	46-271435P1	2	SUPPORT, B R	2	BELT GUIDE, DELRIN
10	2212191	2	CLAMP BELT T	2	DRIVE BELT CLAMP PLATE LPCC
11	46-265089P1	2	SCREW	1	SCREW
12	2212188	2	MOUNTING SUP	2	BEZEL ASSEMBLY SUPPORT LPCC
13	2212187	2	CABLE TRACK	1	CABLE TRACK TO BASEPLATE BKT LPCC
14	2212186	2	CLAMP TAKE-U	1	CABLE TRACK CLAMP LPCC
15	2214666	2	CABLE	1	LCC HD T/R TO HD PREAMP SHIELDED CBL
16	46-264442G3	1	PREAMP BOX	1	1.5T PREAMPLIFIER WITH DIODE FIX
17	46-264988G2	1	CIRCUIT BD	1	LPCC 1.5T HEAD T/R ASSEMBLY
18	2213331-3	2	DOUBLE BEZEL	1	1.5T LPCC BEZEL WITH POKE-A-YOKE
19	2210816	2	CARRIAGE COV	1	LPCC COVER
20	2212193	2	CARRIAGE RAI	1	LEFT RAIL FOR THE COVER ASSEMBLY
21	2212192	2	CARRIAGE RAI	1	RIGHT RAIL FOR THE COVER ASSEMBLY
22	2216528	2	BEZEL BRACKE	1	BEZEL & ELECTRONICS BKT LPCC W/O PA
23	2215931	2	BLANK	1	BLANK FOR DOUBLE BEZEL LPCC_W/O PA
24	2146831	2	WARN LABEL	1	UFI WARNING LABEL - ENGLISH
98	46-208913P36	2	SCREW-INCH	2	SCREW-INCH COUNTERSUNK HEA 1.75 INCH
99	46-208921P36	2	SCREW	2	BRS SLOTTED BND HD SCR 8-32 X 3/8"
102	2223862	2	SPACER	2	BRASS 1/4" TUBULAR SPACER
12923	46-208921P9	2	SCREW	2	BRS SLOTTED BND HD SCR 10-32 X 3/4"
17676	46-208913P7	2	SCREW	4	BRASS SLOTTED FLT HD SCR 8-32 X 1/2"
17716	46-208913P29	2	SCREW	2	BRS SLOTTED FLT HD SCR 10-32 X 1"
17856	46-220184P13	2	SCREW	4	BRS SLOTTED RND HD SCR 6-32 X 3/8"
18062	46-220184P28	2	SCREW	4	BRS SLOTTED RND HD SCR 8-32 X 5/8"
18113	46-220184P41	2	SCREW	11	BRS SLOTTED RND HD SCR 10-32 X 1/2"
18120	46-220184P42	2	SCREW	2	BRS SLOTTED RND HD SCR 10-32 X 5/8"
18129	46-220184P45	2	SCREW	2	BRS SLOTTED RND HD SCR10-32 X 1"
18138	46-220184P47	2	SCREW	4	SCREW__SLOTTED ROUND HEAD BRASS
60436	46-208921P8	2	SCREW	6	BRS SLOTTED BND HD SCR 10-32 X 3/8"
78453	46-221865P3	2	CONNECTOR	1	UG-27B RT-ANGLE N, 50-OHM JACK-PLUG

2-2 1.0T NON-MULTICOIL LPCC ASSEMBLY


Item	Part Number	FRU	Name	Quantity	Description (Remarks)
2	46-282163P2	2	BASE PLATE	1	LPCC CARRIAGE BASEPLATE
3	46-243167P1	2	BOX	12	WHEEL HOUSING, MOLDED NORYL
4	46-243494P1	2	CASTER	6	3/4" ID, 1-3/4 OD, 1/2" WIDE PLASTIC
5	46-243141P1	2	SHLD SCR	4	PLASTIC BEARING MOUNT SHOULDER SCREW
6	46-243178P2	2	BEARING	4	BEARING
7	46-282164P1	2	PLUNGER	1	CARRIAGE RELEASE PLUNGER, DELRIN
8	46-317169P1	2	HOOK	1	CRADLE LATCH HOOK
9	46-271435P1	2	SUPPORT, B R	2	BELT GUIDE, DELRIN
10	2212191	2	CLAMP BELT T	2	DRIVE BELT CLAMP PLATE LPCC
11	46-265089P1	2	SCREW	1	SCREW
12	2212188	2	MOUNTING SUP	2	BEZEL ASSEMBLY SUPPORT LPCC
13	2212187	2	CABLE TRACK	1	CABLE TRACK TO BASEPLATE BKT LPCC
14	2212186	2	CLAMP TAKE-U	1	CABLE TRACK CLAMP LPCC
15	2214666	2	CABLE	1	LCC HD T/R TO HD PREAMP SHIELDED CBL
16	46-321338G3	1	PREAMP BOX	1	1.0T PREAMPLIFIER WITH DIODE FIX
17	46-321336G2	1	CIRCUIT BD	1	LPCC 1.0T HEAD T/R ASSEMBLY
18	2213331-2	2	DOUBLE BEZEL	1	1.0T LPCC DOUBLE BEZEL
19	2210816	2	CARRIAGE COV	1	LPCC COVER
20	2212193	2	CARRIAGE RAI	1	LEFT RAIL FOR THE COVER ASSEMBLY
21	2212192	2	CARRIAGE RAI	1	RIGHT RAIL FOR THE COVER ASSEMBLY
22	2216528	2	BEZEL BRACKE	1	BEZEL & ELECTRONICS BKT LPCC W/O PA
23	2215931	2	BLANK	1	BLANK FOR DOUBLE BEZEL LPCC_W/O PA
98	46-208913P36	2	SCREW-INCH	2	SCREW-INCH COUNTERSUNK HEA 1.75 INCH
99	46-208921P36	2	SCREW	2	BRS SLOTTED BND HD SCR 8-32 X 3/8"
102	2223862	2	SPACER	2	BRASS 1/4" TUBULAR SPACER
12923	46-208921P9	2	SCREW	2	BRS SLOTTED BND HD SCR 10-32 X 3/4"
17676	46-208913P7	2	SCREW	4	BRASS SLOTTED FLT HD SCR 8-32 X 1/2"
17716	46-208913P29	2	SCREW	2	BRS SLOTTED FLT HD SCR 10-32 X 1"
17856	46-220184P13	2	SCREW	4	BRS SLOTTED RND HD SCR 6-32 X 3/8"
18062	46-220184P28	2	SCREW	4	BRS SLOTTED RND HD SCR 8-32 X 5/8"
18113	46-220184P41	2	SCREW	11	BRS SLOTTED RND HD SCR 10-32 X 1/2"
18120	46-220184P42	2	SCREW	2	BRS SLOTTED RND HD SCR 10-32 X 5/8"
18129	46-220184P45	2	SCREW	2	BRS SLOTTED RND HD SCR10-32 X 1"
18138	46-220184P47	2	SCREW	4	SCREW__SLOTTED ROUND HEAD BRASS
60436	46-208921P8	2	SCREW	6	BRS SLOTTED BND HD SCR 10-32 X 3/8"
78453	46-221865P3	2	CONNECTOR	1	UG-27B RT-ANGLE N, 50-OHM JACK-PLUG

This page left blank intentionally.

SECTION 3 – MULTICOIL UPGRADE PROCESS


3-1 MULTICOIL UPGRADE PROCESS

3A – TRANSFER EXISTING PREAMPS AND CONNECT TO BNC INPUT CABLES

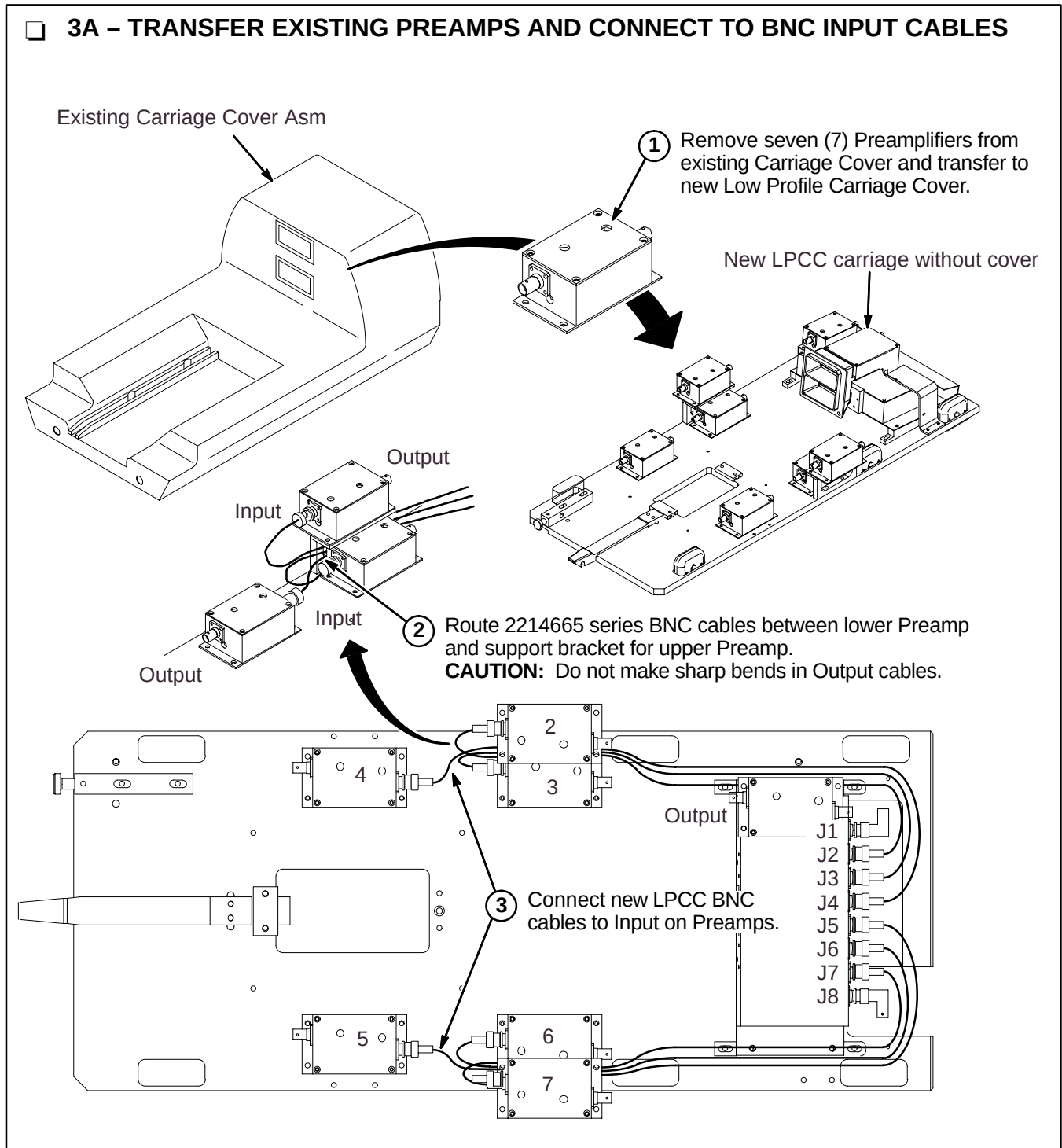

3-1 1.5T MULTICOIL UPGRADE PROCESS

3B – INSTALL AND CONNECT CABLE TRACK BNC OUTPUT CABLES

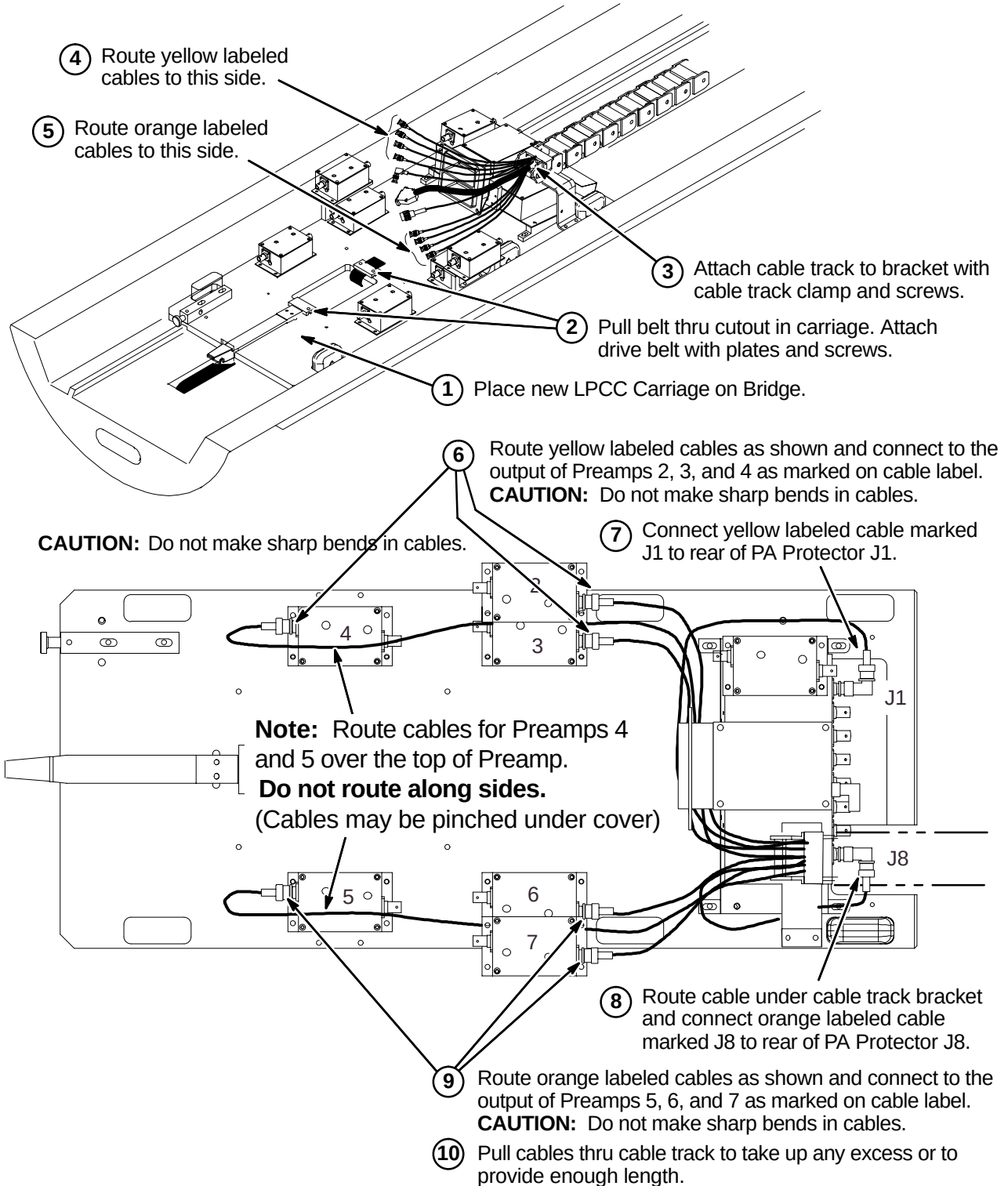


3-1 1.5T MULTICOIL UPGRADE PROCESS

3C – CONNECT RUNS 404, 485 AND 736

① Connect Run 485 marked J4 to Preamp Output.

② Connect Run 736 to J9 on side of Preamp Protection module.
Note: If cable entrance into connector is not as shown, perform steps 4 thru 10.

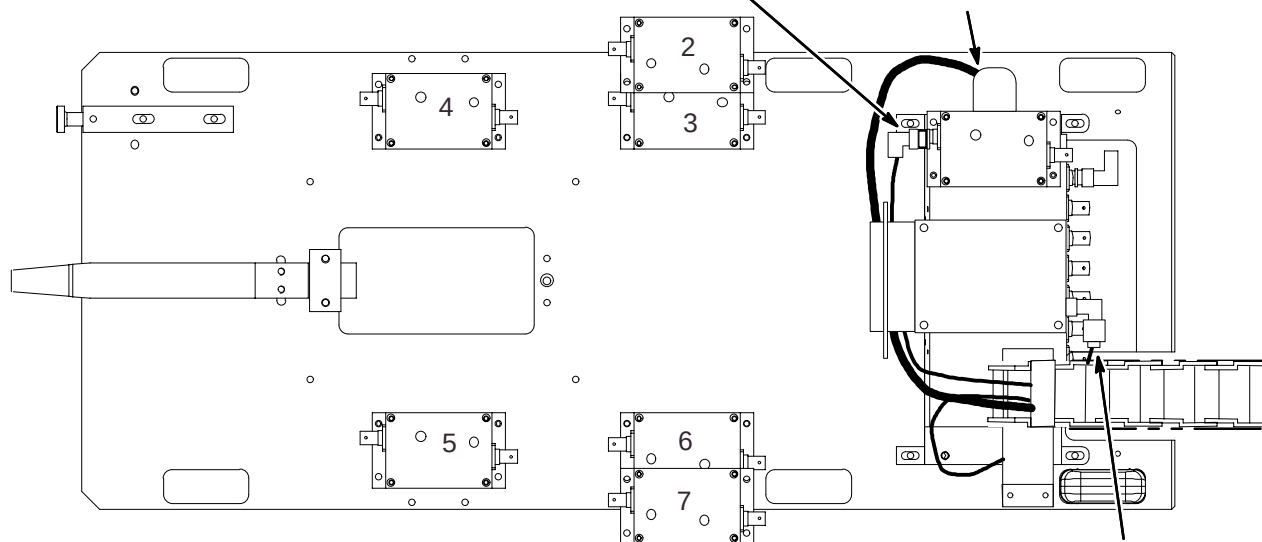

③ Route Run 404 under cable track bracket and connect end marked J6 to Head T/R Assembly.

④ Open connector.

⑤ Remove cover from side hole. Save for installing in top hole.

⑥ Loosen, but do not remove, screw holding cable clamp.

⑦ Remove and save screw.

⑧ Move cable to side entry hole in connector. Rotate cable clamp over cable. Install removed screw from Step 4. Tighten both screws.

⑨ Install cover in top hole.

⑩ Close connector cover.


3-1 1.5T MULTICOIL UPGRADE PROCESS

3D - TY-WRAP CABLES

- ① Route cables under connector bezel. Pull excess cable thru cable track.

- ② Ty-rap cables here.


- ③ Ty-rap cables.

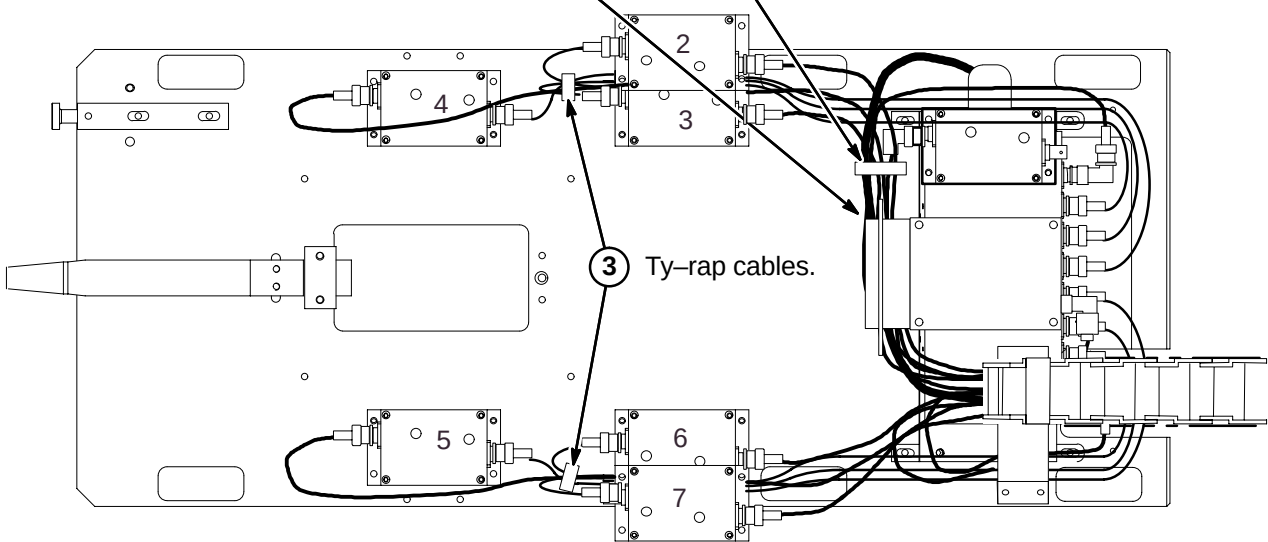

3E - INSTALL CARRIAGE COVER

- ① Install Cover to Carriage Asm.


Note: Slide/position cover forward to clear button and bezel. Then drop down while tucking in insulator. Push cover back to seat on bezel. Then install screws.

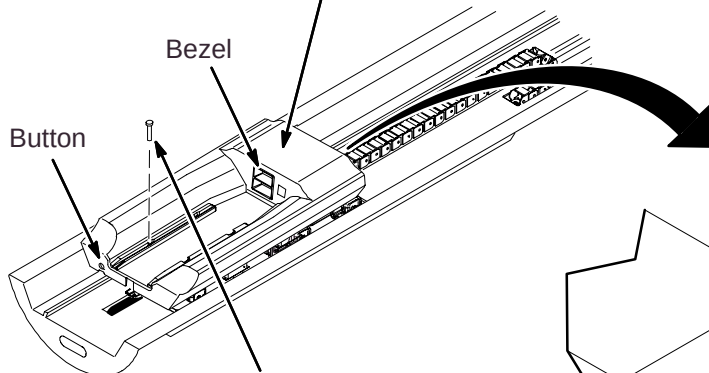

- ② Tuck insulator inside Cover as cover is being installed.


CAUTION: Failure to do this will result in the insulator getting caught in the Cable Track.

- ③ Attach Cover to Carriage using four 10-32 x 1-1/2 Brass round head screws (46-220184P47).

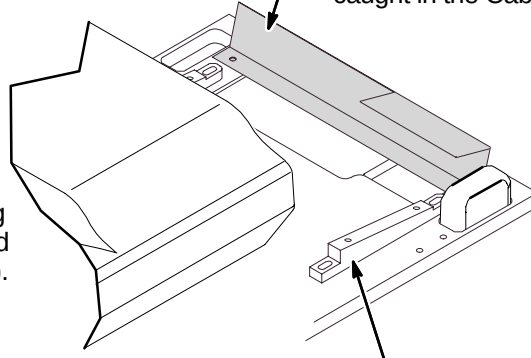

Note: Mounting brackets may be adjusted forward or back for best fit of cover to base.

3-2 1.5T AND 1.0T MULTICOIL LPCC ASSEMBLY

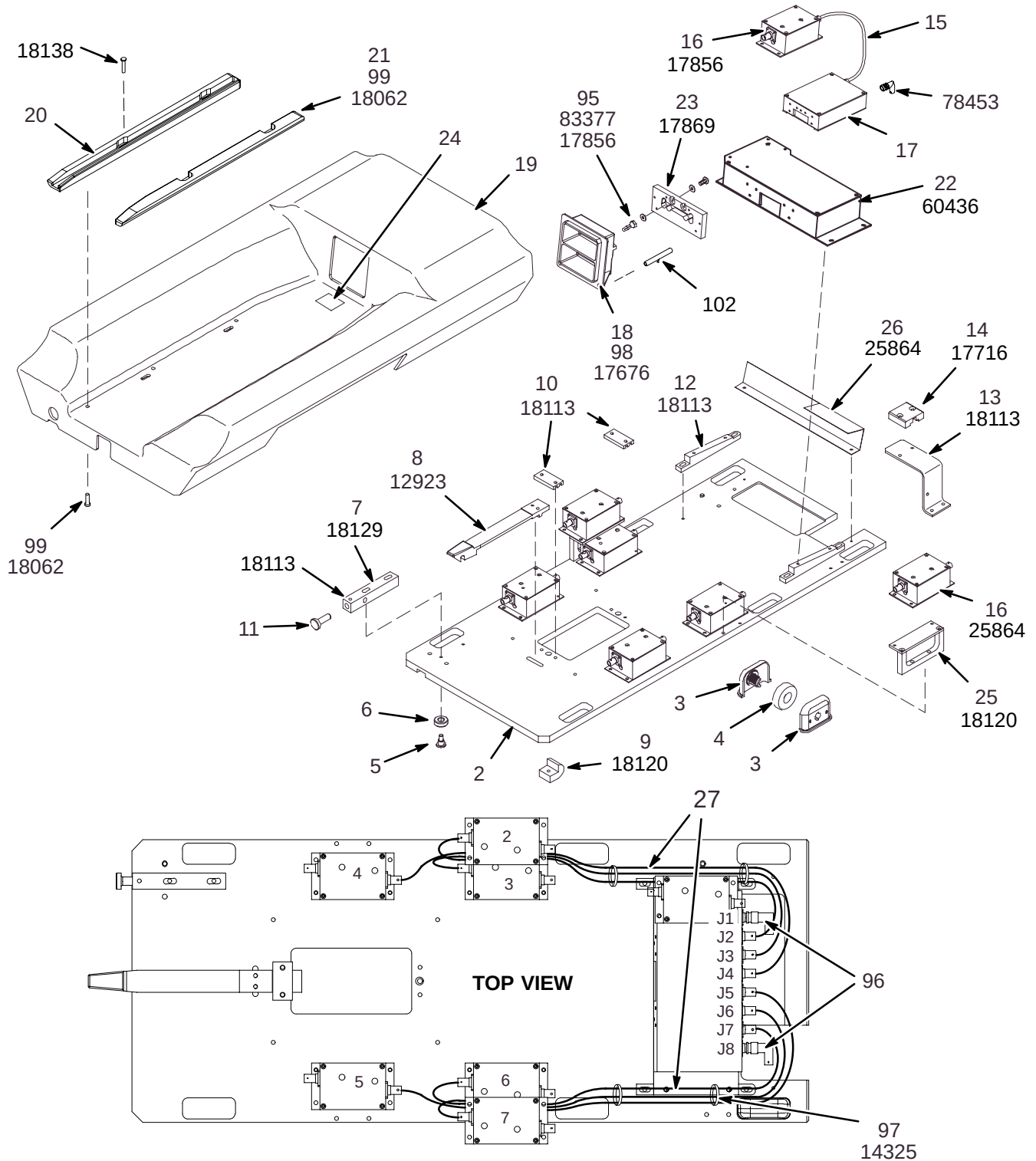

3-2 1.5T MULTICOIL LPCC ASSEMBLY

Item	Part Number	FRU	Name	Quantity	Description (Remarks)
2	46-282163P2	2	BASE PLATE	1	LPCC CARRIAGE BASEPLATE
3	46-243167P1	2	BOX	12	WHEEL HOUSING, MOLDED NORYL
4	46-243494P1	2	CASTER	6	3/4" ID, 1-3/4 OD, 1/2" WIDE PLASTIC
5	46-243141P1	2	SHLD SCR	4	PLASTIC BEARING MOUNT SHOULDER SCREW
6	46-243178P2	2	BEARING	4	BEARING
7	46-282164P1	2	PLUNGER	1	CARRIAGE RELEASE PLUNGER, DELRIN
8	46-317169P1	2	HOOK	1	CRADLE LATCH HOOK
9	46-271435P1	2	SUPPORT, B R	2	BELT GUIDE, DELRIN
10	2212191	2	CLAMP BELT T	2	DRIVE BELT CLAMP PLATE LPCC
11	46-265089P1	2	SCREW	1	SCREW
12	2212188	2	MOUNTING SUP	2	BEZEL ASSEMBLY SUPPORT LPCC
13	2212187	2	CABLE TRACK	1	CABLE TRACK TO BASEPLATE BKT LPCC
14	2212186	2	CLAMP TAKE-U	1	CABLE TRACK CLAMP LPCC
15	2214666	2	CABLE	1	LCC HD T/R TO HD PREAMP SHIELDED CBL
16	46-264442G3	1	PREAMP BOX	7	1.5T PREAMPLIFIER WITH DIODE FIX
17	46-264988G2	1	CIRCUIT BD	1	LPCC 1.5T HEAD T/R ASSEMBLY
18	2213331-3	2	DOUBLE BEZEL	1	1.5T LPCC BEZEL WITH POKE-A-YOKE
19	2210816	2	CARRIAGE COV	1	LPCC COVER
20	2212193	2	CARRIAGE RAI	1	LEFT RAIL FOR THE COVER ASSEMBLY
21	2212192	2	CARRIAGE RAI	1	RIGHT RAIL FOR THE COVER ASSEMBLY
22	46-288510G3	1	CIRCUIT BD	1	LPCC 1.5T PREAMP PROTECTION ASM
23	46-320938P2	2	BEZEL PLATE	1	LPCC 1.5T BEZEL PLATE
24	2146831	2	WARN LABEL	1	UFI WARNING LABEL - ENGLISH
25	2212189	2	PRE AMP MOUN	2	PREAMP SUPPORT ASSEMBLY LPCC W/PA
26	2210732	2	BASE PLATE	1	FISHPAPER INSULATOR LPCC W/PA
27	2214665-9	2	CABLE	1	LPCC PREAMP INPUT CABLE COLLECTOR (SEE PAGE 15)
95	46-282460P1	2	CONNECTOR	2	STD, HEX BRASS UNINSULATED BANANA
96	46-233422P1	2	CONNECTOR	2	BNC JACK-PLUG, RIGHT ANGLE ADAPTER.
97	46-287119P1	2	VIBRATION IS	10	SELF-ADHERE TIE-STRAP ANCHOR, NYLON
98	46-208913P36	2	SCREW-INCH	2	SCREW-INCH COUNTERSUNK HEA 1.75 INCH
99	46-208921P36	2	SCREW	2	BRS SLOTTED BND HD SCR 8-32 X 3/8"
102	2223862	2	SPACER	2	BRASS 1/4" TUBULAR SPACER
12923	46-208921P9	2	SCREW	2	BRS SLOTTED BND HD SCR 10-32 X 3/4"
14325	46-208758P3	2	TY-RAP	12	7.31" X .184" SELF-LOCKING CABLE TIE
17676	46-208913P7	2	SCREW	4	BRASS SLOTTED FLT HD SCR 8-32 X 1/2"
17716	46-208913P29	2	SCREW	2	BRS SLOTTED FLT HD SCR 10-32 X 1"
17856	46-220184P13	2	SCREW	6	BRS SLOTTED RND HD SCR 6-32 X 3/8"
17869	46-220184P15	2	SCREW	2	BRS SLOTTED RND HD SCR 6-32 X 5/8"
18062	46-220184P28	2	SCREW	4	BRS SLOTTED RND HD SCR 8-32 X 5/8"
18113	46-220184P41	2	SCREW	11	BRS SLOTTED RND HD SCR 10-32 X 1/2"
18120	46-220184P42	2	SCREW	6	BRS SLOTTED RND HD SCR 10-32 X 5/8"
18129	46-220184P45	2	SCREW	2	BRS SLOTTED RND HD SCR10-32 X 1"
18138	46-220184P47	2	SCREW	4	SCREW__SLOTTED ROUND HEAD BRASS
25864	46-208921P7	2	SCREW	14	SCREW
60436	46-208921P8	2	SCREW	4	BRS SLOTTED BND HD SCR 10-32 X 3/8"
78453	46-221865P3	2	CONNECTOR	1	UG-27B RT-ANGLE N, 50-OHM JACK-PLUG
83377	46-221449P3	2	WASHER	4	#6 SCREW SIZE (0.138) SILICONE BRONZ

3-2 1.0T MULTICOIL LPCC ASSEMBLY

Item	Part Number	FRU	Name	Quantity	Description (Remarks)
2	46-282163P2	2	BASE PLATE	1	LPCC CARRIAGE BASEPLATE
3	46-243167P1	2	BOX	12	WHEEL HOUSING, MOLDED NORYL
4	46-243494P1	2	CASTER	6	3/4" ID, 1-3/4 OD, 1/2" WIDE PLASTIC
5	46-243141P1	2	SHLD SCR	4	PLASTIC BEARING MOUNT SHOULDER SCREW
6	46-243178P2	2	BEARING	4	BEARING
7	46-282164P1	2	PLUNGER	1	CARRIAGE RELEASE PLUNGER, DELRIN
8	46-317169P1	2	HOOK	1	CRADLE LATCH HOOK
9	46-271435P1	2	SUPPORT, B R	2	BELT GUIDE, DELRIN
10	2212191	2	CLAMP BELT T	2	DRIVE BELT CLAMP PLATE LPCC
11	46-265089P1	2	SCREW	1	SCREW
12	2212188	2	MOUNTING SUP	2	BEZEL ASSEMBLY SUPPORT LPCC
13	2212187	2	CABLE TRACK	1	CABLE TRACK TO BASEPLATE BKT LPCC
14	2212186	2	CLAMP TAKE-U	1	CABLE TRACK CLAMP LPCC
15	2214666	2	CABLE	1	LCC HD T/R TO HD PREAMP SHIELDED CBL
16	46-321338G3	1	PREAMP BOX	7	1.0T PREAMPLIFIER WITH DIODE FIX
17	46-321336G2	1	CIRCUIT BD	1	LPCC 1.0T HEAD T/R ASSEMBLY
18	2213331-2	2	DOUBLE BEZEL	1	1.0T LPCC DOUBLE BEZEL
19	2210816	2	CARRIAGE COV	1	LPCC COVER
20	2212193	2	CARRIAGE RAI	1	LEFT RAIL FOR THE COVER ASSEMBLY
21	2212192	2	CARRIAGE RAI	1	RIGHT RAIL FOR THE COVER ASSEMBLY
22	46-321218G3	1	CIRCUIT BD	1	LPCC 1.0T PREAMP PROTECTION
23	46-320938P3	2	BEZEL PLATE	1	LPCC 1.0T BEZEL PLATE
25	2212189	2	PRE AMP MOUN	2	PREAMP SUPPORT ASSEMBLY LPCC W/PA
26	2210732	2	BASE PLATE	1	FISHPAPER INSULATOR LPCC W/PA
27	2214665-9	2	CABLE	1	LPCC PREAMP INPUT CABLE COLLECTOR (SEE PAGE 15)
95	46-282460P1	2	CONNECTOR	2	STD, HEX BRASS UNINSULATED BANANA
96	46-233422P1	2	CONNECTOR	2	BNC JACK-PLUG, RIGHT ANGLE ADAPTER.
97	46-287119P1	2	VIBRATION IS	10	SELF-ADHERE TIE-STRAP ANCHOR, NYLON
98	46-208913P36	2	SCREW-INCH	2	SCREW-INCH COUNTERSUNK HEA 1.75 INCH
99	46-208921P36	2	SCREW	2	BRS SLOTTED BND HD SCR 8-32 X 3/8"
102	2223862	2	SPACER	2	BRASS 1/4" TUBULAR SPACER
12923	46-208921P9	2	SCREW	2	BRS SLOTTED BND HD SCR 10-32 X 3/4"
14325	46-208758P3	2	TY-RAP	12	7.31" X .184" SELF-LOCKING CABLE TIE
17676	46-208913P7	2	SCREW	4	BRASS SLOTTED FLT HD SCR 8-32 X 1/2"
17716	46-208913P29	2	SCREW	2	BRS SLOTTED FLT HD SCR 10-32 X 1"
17856	46-220184P13	2	SCREW	6	BRS SLOTTED RND HD SCR 6-32 X 3/8"
17869	46-220184P15	2	SCREW	2	BRS SLOTTED RND HD SCR 6-32 X 5/8"
18062	46-220184P28	2	SCREW	4	BRS SLOTTED RND HD SCR 8-32 X 5/8"
18113	46-220184P41	2	SCREW	11	BRS SLOTTED RND HD SCR 10-32 X 1/2"
18120	46-220184P42	2	SCREW	6	BRS SLOTTED RND HD SCR 10-32 X 5/8"
18129	46-220184P45	2	SCREW	2	BRS SLOTTED RND HD SCR10-32 X 1"
18138	46-220184P47	2	SCREW	4	SCREW__SLOTTED ROUND HEAD BRASS
25864	46-208921P7	2	SCREW	14	SCREW
60436	46-208921P8	2	SCREW	4	BRS SLOTTED BND HD SCR 10-32 X 3/8"
78453	46-221865P3	2	CONNECTOR	1	UG-27B RT-ANGLE N, 50-OHM JACK-PLUG
83377	46-221449P3	2	WASHER	4	#6 SCREW SIZE (0.138) SILICONE BRONZ

LPCC PREAMP INPUT CABLES 2214665-9

Item	Part Number	FRU	Name	Quantity	Description (Remarks)
1	2214665-2	2	CABLE	1	LPCC PREAMP #2 INPUT SHIELDED CABLE
2	2214665-3	2	CABLE	1	LPCC PREAMP #3 INPUT SHIELDED CABLE
3	2214665-4	2	CABLE	1	LPCC PREAMP #4 INPUT SHIELDED CABLE
4	2214665-5	2	CABLE	1	LPCC PREAMP #5 INPUT SHIELDED CABLE
5	2214665-6	2	CABLE	1	LPCC PREAMP #6 INPUT SHIELDED CABLE
6	2214665-7	2	CABLE	1	LPCC PREAMP #7 INPUT SHIELDED CABLE