

SECTION 1 – SYSTEM CONFIGURATION

TABLE OF CONTENTS

<u>SECTION</u>	<u>TITLE</u>	<u>PAGE</u>
1-1	INTRODUCTION	1-3
1-2	BASIC SYSTEM	1-4
	1-2-1 System Cooling Equipment Selections	1-11
1-3	SYSTEM OPTIONS	1-14
1-4	FACILITY OPTIONS	1-15
1-5	VAN MANUFACTURERS CONTACT INFORMATION	1-16

1-1 INTRODUCTION

Magnetic Resonance (MR) system Signa® 1.5T TwinSpeed with EXCITE Technology and Signa® 1.5T TwinSpeed are documented throughout this publication.

Note

There are some naming differences between the sales literature documents and the technical documentation. Listed below are the corresponding names. The Pre-Installation document uses the technical document names.

SALES LITERATURE NAME	=	TECHNICAL DOCUMENT NAME
Signa Infinity 1.5T TwinSpeed with EXCITE Technology	=	Signa 1.5T TwinSpeed with EXCITE Technology
Signa Infinity 1.5T TwinSpeed	=	Signa 1.5T TwinSpeed
CXK4 150 Magnet with K4 Technology	=	1.5T LCC (Active Shield) Magnet

Note

Throughout this document references to Signa TwinSpeed apply to both Signa 1.5T TwinSpeed with EXCITE Technology and Signa 1.5T TwinSpeed except where explicitly documented as different.

1-2 BASIC SYSTEM

The basic system for fixed site operation consists of the following major equipment:

- System cooling equipment (See Illustration 1-1 flowchart for system cooling equipment configurations options and relationship to other system equipment/options selections.)


Note

True Mobile system water cooling for the Gradient Coil and the Shield/Cryo Cooler Compressor is provided by the Mobile Van equipment.

- 15 kilogauss (1.5T) LCC (Active Shield) Magnet & Magnet Enclosure
- Magnet and Magnet Accessories
- Shield/Cryo Cooler Compressor Cabinet
- TwinSpeed Gradient (TRM) and RF body coils and 1.5T Quad Head Coil
- Main Disconnect Panel (GE supplied or customer supplied, dependent on Illustration 1-1 flowchart decisions)
- System electronics cabinets: ACGD/PDU Cabinet, 1.5T SRFD II Cabinet, Twin Accessory Cabinet, and the System Cabinet
 - For Signa *TwinSpeed* with EXCITE Technology the System Cabinet contains the Multi Generational Data Acquisition (MGD) chassis, Remote RF chassis, and Driver Module,
 - For Signa *TwinSpeed* the System Cabinet contains the Integrated System Electronics (ISE) chassis and Transient Noise Suppression Module
- Operator Workspace (OW) equipment
- Patient Transport Table
- Patient accessories

The tables in this section list the catalogs for the hardware equipment which comprise the Signa 1.5T *TwinSpeed* system. Illustrations 1-2 and 1-3 show the major equipment of the system. Refer to Section **1-2-1 System Cooling Equipment Selections** for explanation of available configurations.

1-2 BASIC SYSTEM (Continued)

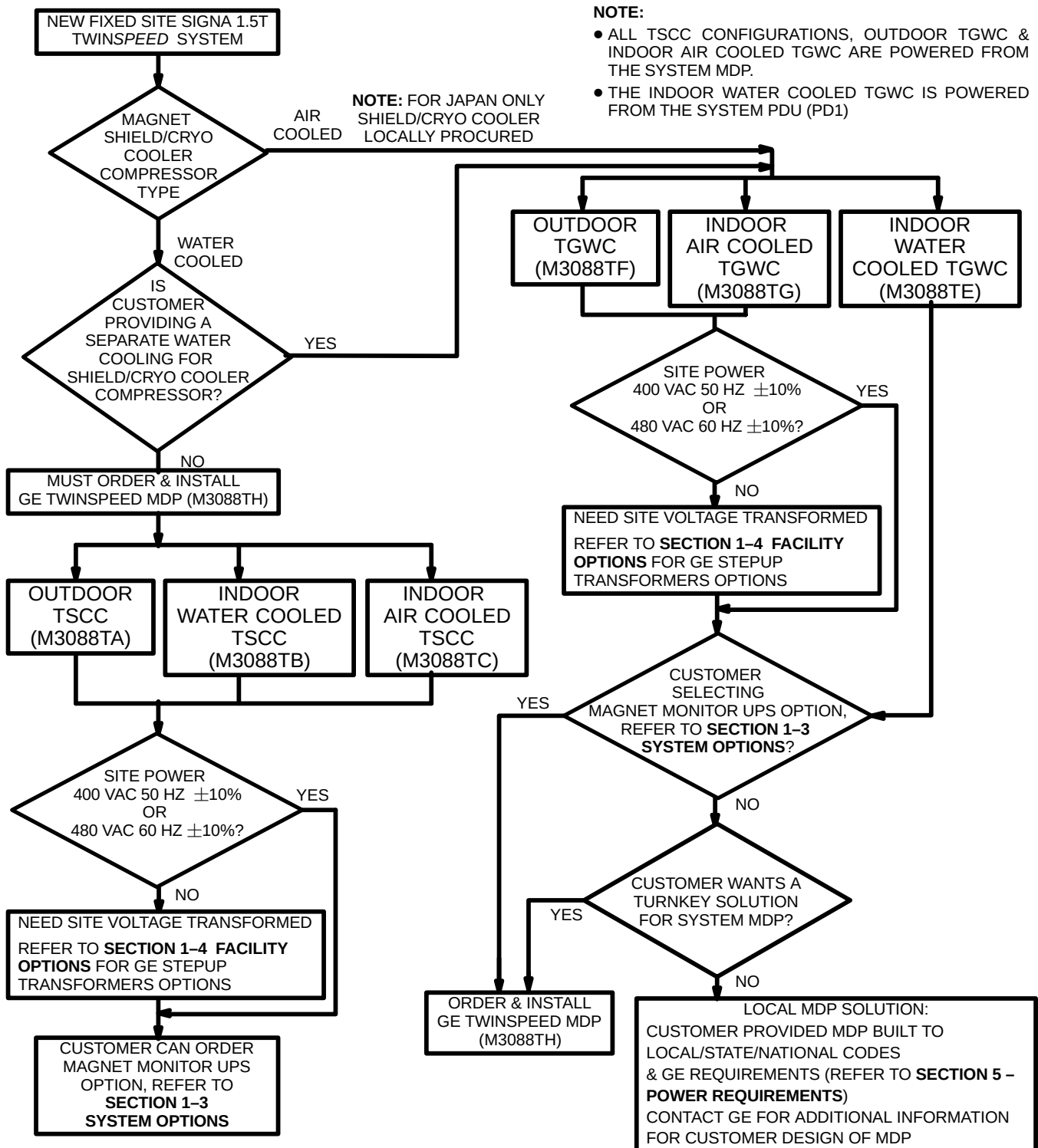

FIXED SITE SYSTEM COOLING EQUIPMENT & MAIN DISCONNECT PANEL SELECTIONS

ILLUSTRATION 1-1

1-2 BASIC SYSTEM (Continued)

TABLE 1-1
SYSTEM ELECTRONICS MAJOR EQUIPMENT

CATALOG	DESCRIPTION
M3000TC Signa 1.5T TwinSpeed with EXCITE Technology	• System Cabinet consisting of MGD chassis, Remote RF chassis, and Driver Module. • ACGD/PDU Cabinet containing ACGD Gradient drivers and Power Distribution Unit module with unregulated transformer 200/380/400/415/480 Volt, 50/60 Hz with power filter. • 1.5T SRFD II Cabinet containing the RF Amplifiers, RF Power Monitor, power supplies for the Magnet Enclosure system components. • Twin Accessory Cabinet containing the Gradient Switch and the vacuum pump for Quiet Technology Magnet Enclosure. • Operator Workspace (OW) equipment: Octane 2 Computer, Workspace Cabinet, Mouse and Mouse Pad, LCD panel, and chair. Note, refer to Tables 1-6, 1-7, and 1-8 for catalog choices required to complete Operator Workspace equipment. • Pneumatic Patient Alert System. • 1.5T Quad Head Coil • Patient Transport Table and cradle. • Patient accessories such as: a phantom kit, patient log book, head cushion and sponges, chin and forehead straps, body wedges, knee cushions, and security/restraint straps. • Gating accessories which include: patient cardiac leads, peripheral gating probe, and respiratory bellows.
M3000TA Signa 1.5T TwinSpeed	• System Cabinet consisting of Integrated Systems Electronics subsystem, Universal Combined Exciter/Receiver (UCERD), Phase Array hardware integrated, and with Fast Receiver. • ACGD/PDU Cabinet containing ACGD Gradient drivers and Power Distribution Unit module with unregulated transformer 200/380/400/415/480 Volt, 50/60 Hz with power filter. • 1.5T SRFD II Cabinet containing the RF Amplifiers, RF Power Monitor, power supplies for the Magnet Enclosure system components. • Twin Accessory Cabinet containing the Gradient Switch and the vacuum pump for Quiet Technology Magnet Enclosure. • Operator Workspace (OW) equipment: Octane Computer, Workspace Cabinet, Mouse and Mouse Pad, LCD panel, and chair. Note, refer to Tables 1-6, 1-7, and 1-8 for catalog choices required to complete Operator Workspace equipment. • Pneumatic Patient Alert System. • 1.5T Quad Head Coil • Patient Transport Table and cradle. • Patient accessories such as: a phantom kit, patient log book, head cushion and sponges, chin and forehead straps, body wedges, knee cushions, and security/restraint straps. • Gating accessories which include: patient cardiac leads, peripheral gating probe, and respiratory bellows.

1-2 BASIC SYSTEM (Continued)

TABLE 1-2
1.5T MAGNET SUBSYSTEM (SELECT ONE)

CATALOG	DESCRIPTION
M3060TA	15 kilogauss (1.5T) LCC (Active Shield) Magnet for Fixed Site or Transportable/Relocatable Mobile (See Note 1). Includes Quiet Technology Enclosure partially installed on magnet, Magnet Rundown Unit, Magnet Monitor (performs cryogen meter function and magnet pressure control), and Global Modem. (See Notes 2 & 3)
M3061TA	15 kilogauss (1.5T) LCC (Active Shield) Magnet for True Mobile (See Note 1). Includes Quiet Technology Enclosure partially installed on magnet, Magnet Rundown Unit, Magnet Monitor (performs cryogen meter function and magnet pressure control), UPS for Magnet Monitor equipment, and Global Modem . (See Note 2).
Note 1 For Mobiles (Transportable, Relocatable, and True Mobile) contact your selected Van Manufacturer for site pre-installation, refer to Section 1-5, VAN MANUFACTURERS CONTACT INFORMATION. 2 M1060JW Shield/Cryo Cooler Compressor (water cooled) must be ordered for LCC Magnet. The Shield/Cryo Cooler Compressor is not part of the LCC Magnet Catalogs but is required at magnet delivery. Refer to Section 5 POWER REQUIREMENTS and Section 4-4 Water Cooling Requirements for details of compressor 24 hour power and water cooling needs. 3 If M1060TA Magnet is to be shipped by air then must order M1060SC CX Magnet Air Transport Crate (not required within Continental USA).	

TABLE 1-3
LCC MAGNET 2.5 METER CEILING HEIGHT KIT (OPTIONAL)

CATALOG	DESCRIPTION
M1060SR	LCC Magnet Low Ceiling Height (2.5 Meter) Siting Option Note: This option reduces the minimum required ceiling height for the Magnet Room to 2.5 meters from 2.67 meters. Check with GE Installation Specialist to determine if this is a required option.

TABLE 1-4
INTEGRATED BODY MODULE (SELECT ONE)

CATALOG	DESCRIPTION
M3085TA	1.5T Advanced 60 cm Twin Integrated Body Module (TRM)

TABLE 1-5
MAIN DISCONNECT PANEL (MDP) (SEE TABLE NOTE)

CATALOG	DESCRIPTION
M3088TH	Signa TwinSpeed Main Disconnect Panel
Note 1 See Illustration 1-1 flowchart to determine when GE MDP is to be utilized or customer supplied MDP which meets requirements can be considered.	

TABLE 1-6
OPERATOR WORKSPACE MONITOR (SELECT ONE)

CATALOG	DESCRIPTION
M1000NT	18 inch LCD Color Monitor

1-2 BASIC SYSTEM (Continued)

TABLE 1-7
OPERATOR WORKSPACE TABLE (ONE REQUIRED)

CATALOG	DESCRIPTION
M1000MW	Table for Operator Workspace See Note 1
Note 1 OW Table is an integral part of the regulatory approved system. OW Table provides mounting for several assemblies (e.g. OW Interface Module, OW Power Distribution Box, Modem, DASM) and cable routing for OW interconnects.	

TABLE 1-8
COUNTRY KITS (SELECT ONE)

CATALOG	DESCRIPTION
M1000MN	English Keyboard
M1000MP	French Keyboard
M1000MR	German Keyboard
M1000MS	Scandinavian Keyboard
M1000NH	Italian Keyboard
M1000NJ	Portuguese Keyboard
M1000NK	Spanish Keyboard

TABLE 1-9
RECON PROCESSOR WITH BAM MEMORY (SELECT ONE)

CATALOG	DESCRIPTION
FOR SIGNA 1.5T TWINSPEED WITH EXCITE TECHNOLOGY	
M3000RD	EXCITE ReFlex100
M3000RE	EXCITE ReFlex200
M3000RF	EXCITE ReFlex400
FOR SIGNA 1.5T TWINSPEED	
M3000RA	ReFlex100 – Hex Power PC Processor with four 128 MB Bulk Acquisition Memory (BAM) boards (512 MB)

TABLE 1-10
FILMING INTERFACE (SEE NOTE 1 IN TABLE)


CATALOG	DESCRIPTION
M1000MJ	Analog DASM
M1000MK	Digital DASM
Note 1 Must choose one Filming DASM Board unless DICOM Print will be used exclusively for software filming to DICOM print peripheral devices.	

1–2 BASIC SYSTEM (Continued)

TABLE 1–11
SITE COLLECTOR KIT (SELECT ONE)

CATALOG	DESCRIPTION
FOR SIGNA 1.5T TWINSPEED WITH EXCITE TECHNOLOGY	
M3043TC Refer to Notes 1 & 2	● LCC (Active Shield) Magnet Fixed Site Quiet Technology Enclosure ● Penetration Panel with Phase Array hardware integrated and Penetration Panel Covers ● Fixed Site system interconnect cables including Phased Array Option ● Fixed Site Power Distribution Unit cables ● SPT Phantom Set with Shipping/Storage Cart
M3053TC	● LCC (Active Shield) Magnet True Mobile Site Quiet Technology Enclosure ● Penetration Panel with Phase Array hardware integrated ● True Mobile Site system interconnect cables including Phased Array Option ● True Mobile Site Power Distribution Unit cables ● SPT Phantom Set
M3063TC	● LCC (Active Shield) Magnet Transportable/Relocatable Mobile Site Quiet Technology Enclosure ● Penetration Panel with Phase Array hardware integrated ● Transportable/Relocatable Mobile Site system interconnect cables including Phased Array Option ● Transportable/Relocatable Mobile Site Power Distribution Unit cables ● SPT Phantom Set with Shipping/Storage Cart
FOR SIGNA 1.5T TWINSPEED	
M3043TA Refer to Notes 3 & 4	● LCC (Active Shield) Magnet Fixed Site Quiet Technology Enclosure ● Penetration Panel with Phase Array hardware integrated and Penetration Panel Covers ● Fixed Site system interconnect cables including Phased Array Option ● Fixed Site Power Distribution Unit cables ● SPT Phantom Set with Shipping/Storage Cart
M3053TA	● LCC (Active Shield) Magnet True Mobile Site Quiet Technology Enclosure ● Penetration Panel with Phase Array hardware integrated ● True Mobile Site system interconnect cables including Phased Array Option ● True Mobile Site Power Distribution Unit cables ● SPT Phantom Set
M3063TA	● LCC (Active Shield) Magnet Transportable/Relocatable Mobile Site Quiet Technology Enclosure ● Penetration Panel with Phase Array hardware integrated ● Transportable/Relocatable Mobile Site system interconnect cables including Phased Array Option ● Transportable/Relocatable Mobile Site Power Distribution Unit cables ● SPT Phantom Set with Shipping/Storage Cart
Note 1 EXCITE Technology Operator Workspace Long Cables (M1090ZP) is available to order only if the GE Installation Specialist determines the standard length cables (supplied in listed catalogs) for interconnects from OW to other system components will not work at the site. M1090ZP provides long length interconnects for OW to PDU, OW to System Cabinet, and OW to Penetration Panel ONLY.	
2 Signa TwinSpeed with EXCITE Technology Site Collector Kit with Long Cables (M3044TC) is available to order only if the GE Installation Specialist determines the standard length cables in M3043TC will not work at the site. M3044TC includes OW Long Cables with the long cables for the Equipment Room.	
3 Operator Workspace Long Cables (M1090ZN) is available to order only if the GE Installation Specialist determines the standard length cables (supplied in listed catalogs) for interconnects from OW to other system components will not work at the site. M1090ZN provides long length interconnects for OW to PDU, OW to System Cabinet, and OW to Penetration Panel ONLY.	
4 Signa TwinSpeed Site Collector Kit with Long Cables (M3044TA) is available to order only if the GE Installation Specialist determines the standard length cables in M3043TA will not work at the site. M3044TA includes OW Long Cables with the long cables for the Equipment Room.	

1-2 Basic System (Continued)

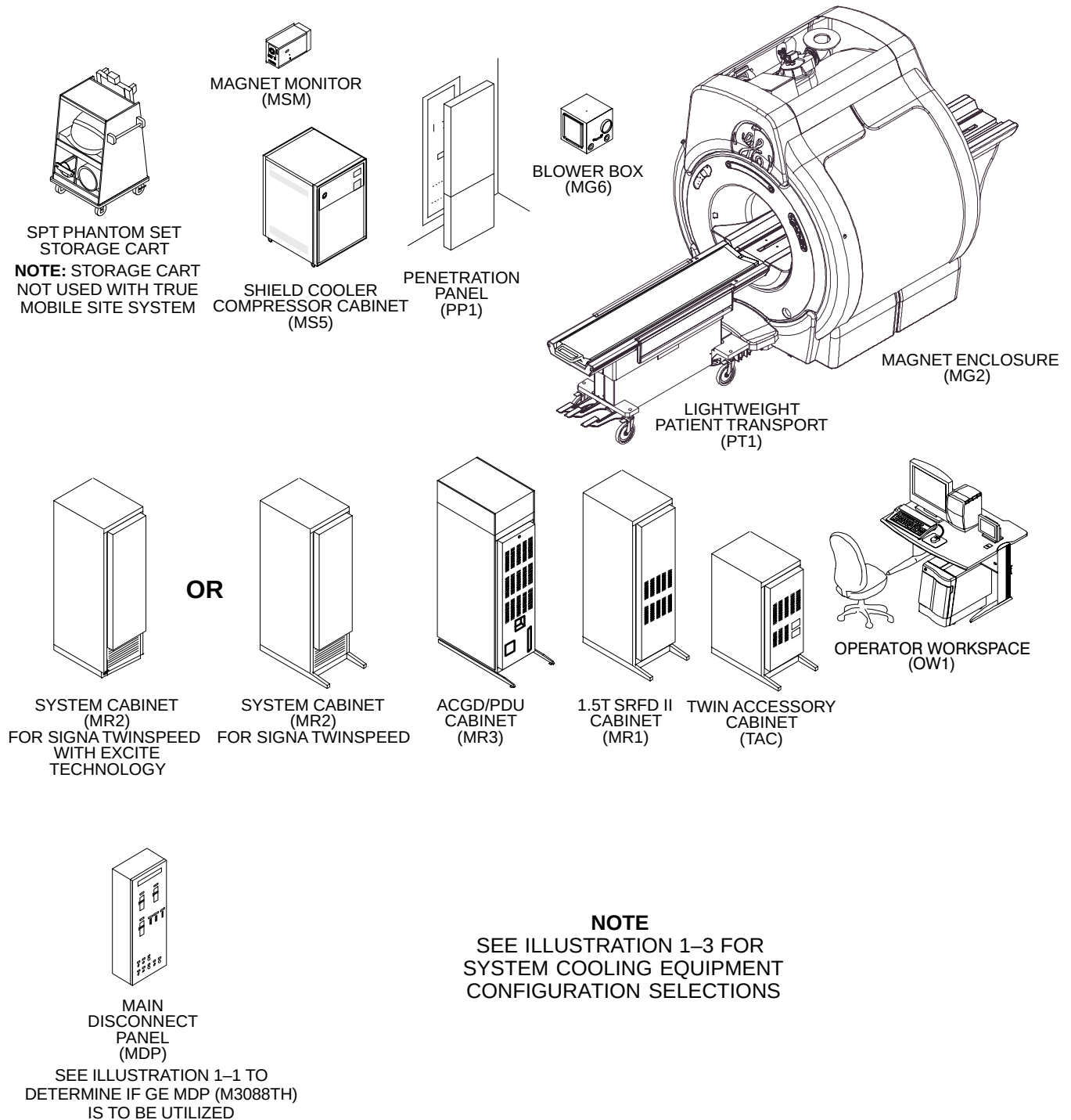

SIGNA 1.5T TWINSPEED SYSTEM
ILLUSTRATION 1-2

1–2–1 System Cooling Equipment Selections

The Signa 1.5T TwinSpeed with EXCITE Technology and the Signa 1.5T TwinSpeed system requires water cooling for the LCC Magnet Shield/Cryo Cooler Compressor Cabinet and the TRM Gradient Coil located in the Magnet bore. The MR Engineering Team together with Ellis & Watts International Inc. have designed several different water cooling cabinet options to optimize siting power requirements and monthly maintenance costs for the customer. The configuration options include:

- **Twin System Cooling Cabinet (TSCC)**
TSCC units provide water cooling for the LCC Magnet Shield/Cryo Cooler Compressor Cabinet and the Gradient Coil located in the Magnet bore. See Table 1–12 for listing of TSCC offerings.
- **Twin Gradient Water Cooler (TGWC)**
TGWC units provide water cooling for ONLY the Gradient Coil located in the Magnet bore. **When a TGWC is selected then the customer site must still provide water cooling for the LCC Magnet Shield/Cryo Cooler Compressor Cabinet.** See Table 1–13 for listing of TGWC offerings.

Note

See Illustration 1–1 flowchart for system cooling equipment configurations options and relationship to other system equipment/options selections.

TSCC Configurations

Refer to Table 1–12 for list of TSCC catalogs offerings.

TABLE 1–12
TWIN SYSTEM COOLING CABINET (TSCC) SEE TABLE NOTE 1

CATALOG	DESCRIPTION
M3088TA	• Outdoor Air Cooled TSCC (Water-to-Air Chiller) 50/60 Hz for Gradient Coil Cooling & Shield/Cryo Cooler Compressor cooling See Note 2 <i>TSCC installed external to building with Remote Control Panel installed in Equipment Room. The most cost effective for the customer.</i>
M3088TB	• Indoor Water Cooled TSCC (Water-to-Water Chiller) 50/60 Hz for Gradient Coil Cooling & Shield/Cryo Cooler Compressor cooling See Note 3 <i>Installed in the Equipment Room, low additional Air Conditioning requirements.</i>
M3088TC	• Indoor Air Cooled TSCC (Water-to-Air Chiller) 50/60 Hz for Gradient Coil Cooling & Shield/Cryo Cooler Compressor cooling <i>Installed in the Equipment Room, increases Air Conditioning requirements.</i>
Note 1 The TSCC requires utilization of the GE MDP per Illustration 1–1. Therefore must order and install M3088TH TwinSpeed Main Disconnect Panel. 2 M3088TA is standard cooling equipment configuration of Transportable/Relocatable Mobile systems. 3 The water cooled configuration System Cooling Cabinets require customer supplied water for cooling. Refer to Section 4–4–3 Water Cooled System Cooling Cabinet Configuration Requirements for water cooling specifications.	

Note

True Mobile system water cooling for the Gradient Coil and the Shield/Cryo Cooler Compressor is provided by the Mobile Van equipment.

1-2-1 System Cooling Equipment Selections (Continued)**TGWC Configurations**

Refer to Table 1-13 for list of TGWC offerings.

TABLE 1-13
TWIN GRADIENT WATER COOLER (TGWC) See Note 1


CATALOG	DESCRIPTION
M3088TF	Outdoor Air Cooled TGWC (Water-to-Air Chiller) 50/60 for Gradient Coil Cooling <i>TGWC installed external to building with Remote Control Panel installed in Equipment Room. The most cost effective for the customer with separate cooling available for the Shield/Cryo Cooler Compressor.</i>
M3088TE	Indoor Water Cooled TGWC (Water-to-Water Heat Exchanger) 50/60 Hz for Gradient Coil Cooling See Note 2 <i>Unit and Remote Control Panel installed in the Equipment Room, low additional Air Conditioning requirements.</i>
M3088TG	Indoor Air Cooled TGWC (Water-to-Air Chiller) 50/60 Hz for Gradient Coil Cooling <i>Installed in the Equipment Room, increases Air Conditioning requirements.</i>
Note 1 When a TGWC is selected then the customer site must still provide water cooling for the LCC Magnet Shield/Cryo Cooler Compressor Cabinet. 2 The water cooled configuration of TGWC requires customer supplied water for cooling. Refer to Section 4-4-3 Water Cooled System Cooling Cabinet Configuration Requirements for water cooling specifications.	

1-2-1 System Cooling Equipment Selections (Continued)


TWIN SYSTEM COOLING CABINET (TSCC)

OUTDOOR TSCC
(M3088TA)


REMOTE
CONTROL
PANEL
(RCP)

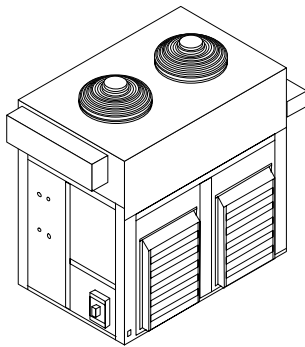

INDOOR WATER
COOLED TSCC
(M3088TB)

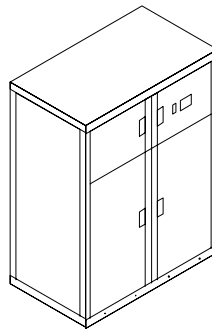

INDOOR AIR
COOLED TSCC
(M3088TC)

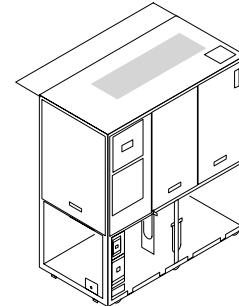


TWIN GRADIENT WATER COOLER
(TGWC)

CAUTION!!

WHEN A TGWC IS SELECTED THEN THE CUSTOMER SITE MUST PROVIDE WATER COOLING FOR THE LCC MAGNET SHIELD/CRYO COOLER COMPRESSOR CABINET.


OUTDOOR UNIT
(M3088TF)

REMOTE
CONTROL
PANEL
(RCP)

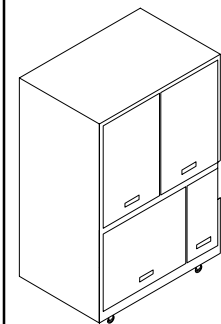


INDOOR WATER
COOLED UNIT
(M3088TE)

REMOTE
CONTROL
PANEL
(RCP)

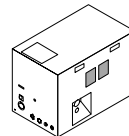

INDOOR AIR
COOLED UNIT
(M3088TG)

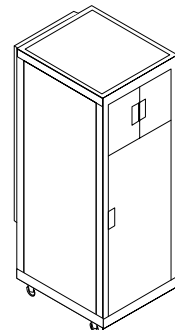

SYSTEM COOLING EQUIPMENT SELECTIONS
ILLUSTRATION 1-3

1-3 SYSTEM OPTIONS

This section lists options for the Signa 1.5T TwinSpeed system which have site preparation impact.

- **For Signa 1.5T TwinSpeed with EXCITE Technology systems ONLY:**

- EXCITE 8 Channel Option (M3033CA)

- 700 Volt Amp UPS for MR Magnet Monitor Equipment (E4504AG)

Note

The GE MDP (M3088TH) must be ordered and installed when the UPS for MR Magnet Monitor (E4504AG) is installed. The GE MDP provides Emergency Off control of the output of the UPS for system safety.

Note

The UPS for MR Magnet Monitor is standard equipment included in the Magnet for True Mobile (M3061TA).

- 2nd Order Resistive Shim Option for Signa 1.5T TwinSpeed (M3080TA). The 2nd Order Shim Option correct inhomogeneities induced by the patient and restores the original homogeneity of the magnet. The Resistive Shim Module is installed in the Twin Accessory Cabinet.
- Advantage Workstation (AW). Refer to Direction 2261309-100 *Advantage Workstation Ultra60 Pre-Installation*.

Note

LCD monitor is recommended for the AW to be located near the Signa system Operator Workspace due to gauss field proximity limits.

- VCR Interface Kit for Octane Computer (M1090TZ)
- Bar Code Reader (M1090PM) for use in conjunction with HIS/RIS software option.
- Additional Patient Transport table with cradle for WideOpen Enclosure (M1090TL).
- Multi-Nuclear Spectroscopy (MNS) option (M1033MD) provides MNS software and hardware. The hardware includes MNS Cabinet containing RF Broadband Amplifier, system interconnect cables, modification hardware for the basic 1.5T system, P31 Flexible Coil, and MR Spectroscopy Head Phantom.
- Advantage Windows Workstation
- **For Signa 1.5T TwinSpeed with EXCITE Technology systems ONLY:** BrainWave Option (M1033BM)
The BrainWave option provides an integrated set of software and hardware stimulus components for the Signa/AW RTIP application to produce brain activation images from MR Blood Oxygen Level Dependent (BOLD) scan data. Stimulus devices include audio (headphones), visual (glasses), and push-button response boxes held by the patient. Stimulus software resides on the AW (AW4.0 or later required) used by RTIP and on a new PC in the BrainWave Cabinet located in the Equipment Room.
- Other hard copy devices and patient accessories.
- Cryogen refill service.
- Various GE Medical Systems Group service contracts.

1-4 FACILITY OPTIONS

- The GE MDP, all configurations of TSCC, the Outdoor TGWC, and the Indoor Air Cooled TGWC require high voltage power. If system input voltage configuration does not meet specifications in **SECTION 5 – POWER REQUIREMENTS** then customer provided transformer may be needed. Refer to Table 1–14 for low voltage step up transformer options.

TABLE 1–14
LOW VOLTAGE (200 or 208 VOLTS) STEP UP TRANSFORMER OPTIONS

CATALOG See Note 1	DESCRIPTION
E4500AS	• 150 KVA 208 to 480Y277 Volt, 60 Hz transformer
R4500BD	• 150 KVA 200 to 400Y230 Volt, 50/60 Hz transformer
Note 1 A step-up transformer is required if the site input voltage configuration does not meet specifications in SECTION 5 – POWER REQUIREMENTS . The GE Main Disconnect Panel (MDP) & Twin System Cooling Cabinet (TSCC) or Twin Gradient Water Cooler (TGWC) configurations require high voltage power.	

- Direct current (DC) lighting controller for the magnet room:
 - E4503AD 20 Amp Maximum Constant Lighting Level System, surface/semi-flush mount
 - E4503AF 20 Amp Maximum Variable Lighting Level System, surface/semi-flush mount
 - E4503AW 28 Amp Maximum Constant Lighting Level System, surface/semi-flush mount
 - E4503AY 28 Amp Maximum Variable Lighting Level System, surface/semi-flush mount.

Note

DC Lighting Controller 50 Hz designs are available by special order, for more information contact:

Dwight Gilbert, GE Supply 414–527–6638 Dwight.Gilbert@supply.ge.com

Lou Hernandez, GE Industrial Systems 262–797–4910 lou.hernandez.indsys.ge.com

- Signature 5000 Series 3 UPS 100KVA (E4502FB) provides reliable, clean, constant voltage power for the complete MR system. The use of a full system UPS enables the system imaging to be completed after loss of supply power and allows for saving of valuable data and orderly system shutdown. Also recommend installing a 100KVA UPS Bypass Panel (E4504CG) which feeds power to the UPS in the normal mode and enables the imaging system to operate when the UPS is in manual bypass mode for routine servicing of the UPS.

- Signa System Seismic Anchorage Service (R4390JA).

Note

Magnet Seismic anchoring is the customer's responsibility to coordinate magnet mounting methods with the RF shielded room vendor to prevent RF leaks and secondary grounding problems. Refer to Section 7–7 ANCHOR HARDWARE REQUIREMENTS FOR MR EQUIPMENT INSIDE RF SHIELDED ROOM for details.

- Oxygen Monitor Kit (M1060KM) which includes Oxygen Monitor and Remote Oxygen Monitor Module.

Note

The Oxygen Monitor does not bear a CE monogram and therefore may not be acceptable in all European countries.

1-5 VAN MANUFACTURERS CONTACT INFORMATION

Pre-Installation requirements may vary for specific Mobile configuration (True Mobile, Transportable, or Relocatable) and between Van Manufacturers. Listed below is contact information for GE Van Manufacturers.

- Ellis & Watts
4400 Glen Willow Lake Lane
Batavia, OH 45103
USA
Telephone: 513-752-9000
FAX: 513-752-4983
- PDC Facilities
700 Walnut Ridge Drive, PO Box 900
Hartland, WI 53029-0900
USA
Telephone: 262-367-7700
FAX: 262-367-7744
- SMIT Mobile Equipment
Buys Ballotstraat 6
3261 La Oud-Beijerland
Holland
Telephone: (31) 186-6-14322
FAX: (31) 186-6-19367
- AK Specialty Vehicles
16745 South Lathrop Avenue
Harvey, IL 60426
USA
Telephone: 708-596-5066
FAX: 708-596-2480
AK Specialty Vehicles is also the contact for used Calumet Coach trailers.