

TABLE OF CONTENTS

TABLE OF CONTENTS 1

1- REAR END BELL REMOVAL 2

2- REAR END BELL INSTALLATION..... 7

REVISION HISTORY 8

Description - Procedures for removing the rear end bell.

Note

If you are replacing the Rear End Bell the gradient cables must be cut. They must be re-crimped and re-terminated.

Parts needed are:

Compression splicer, part number **2218149**. Quantity needed is 12

Cold Shrink, part number **2218148**. Quantity needed is 12.

Tool needed:

Gradient cable crimp tool part number **2135776**

Heavy duty cutters for gradient cables

1- REAR END BELL REMOVAL


1. Remove the front section of the bridge. Rear section can remain attached on the rear pedestal. See procedure for Bridge Removal Twin Enclosure (ME2REA13.DOC).
2. Lock out and tag out circuit breaker MR1 A14 CB6 at rear of RF/Pen Cabinet (MR1) on magnet enclosure power supply (MEPS) module.

For **RF/Pen II** and **RF/PDU Cabinet**: Lock out and tag out circuit breaker MR1 A20 CB4 at rear of RF/Pen II and RF/PDU Cabinet on System Support Module (SSM).

(Refer to *Procedure For Safety: Section 6.*)


3. Verify that voltage between MR1 A14 TP4 (+) and TP13 (–) is 0 Vdc.
4. For **RF/Pen II** and **RF/PDU Cabinet**, verify that voltage between MR1 A20 TP4 (+) (SRI+24V) and TP13 (–) (SRI RTN) is 0 Vdc.

5. Disconnect the microphone cable from the intercom board located on rear pedestal. See Illustration 1-1. The comfort module must be removed to replace the rear end bell.

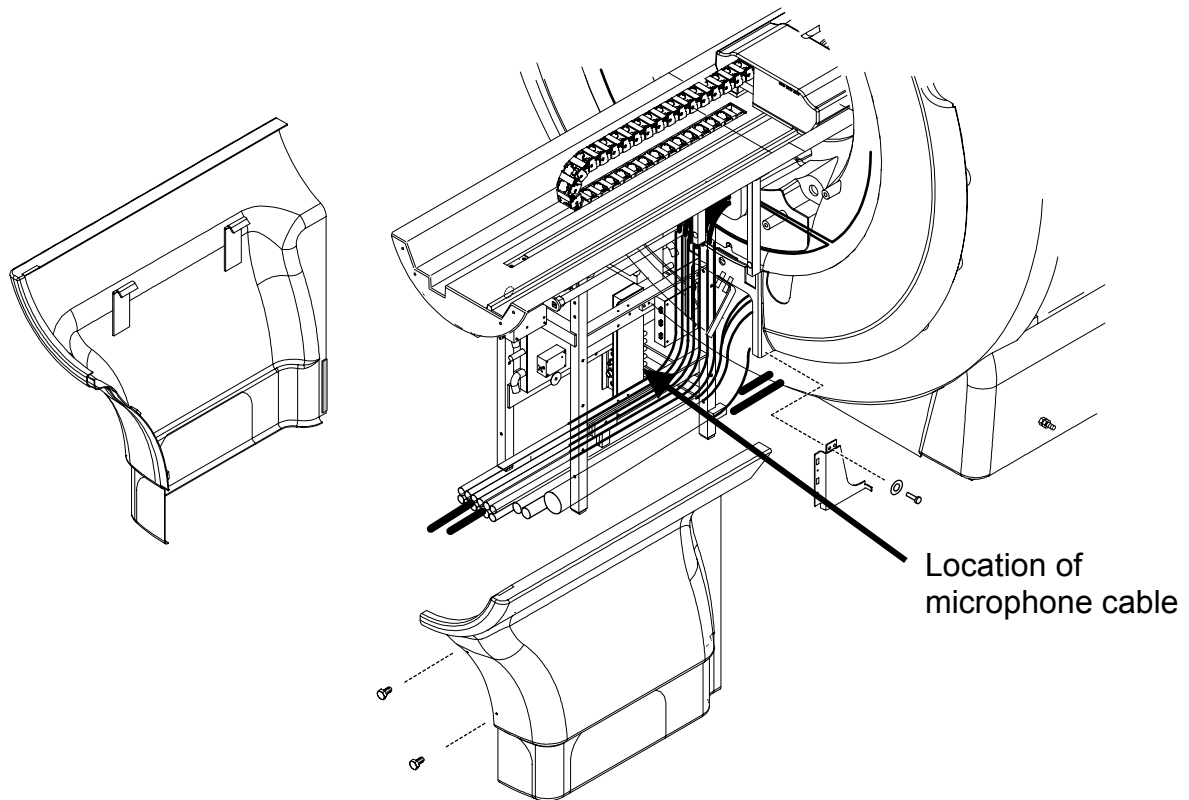


REAR COVERS
ILLUSTRATION 1-1

6. Unclip the light source at the lower right section of the magnet. See Illustration 1-2.
7. Route the cable that runs from the light source up to the comfort module, as this cable will remain connected to the patient lights during the rear end bell replacement process.

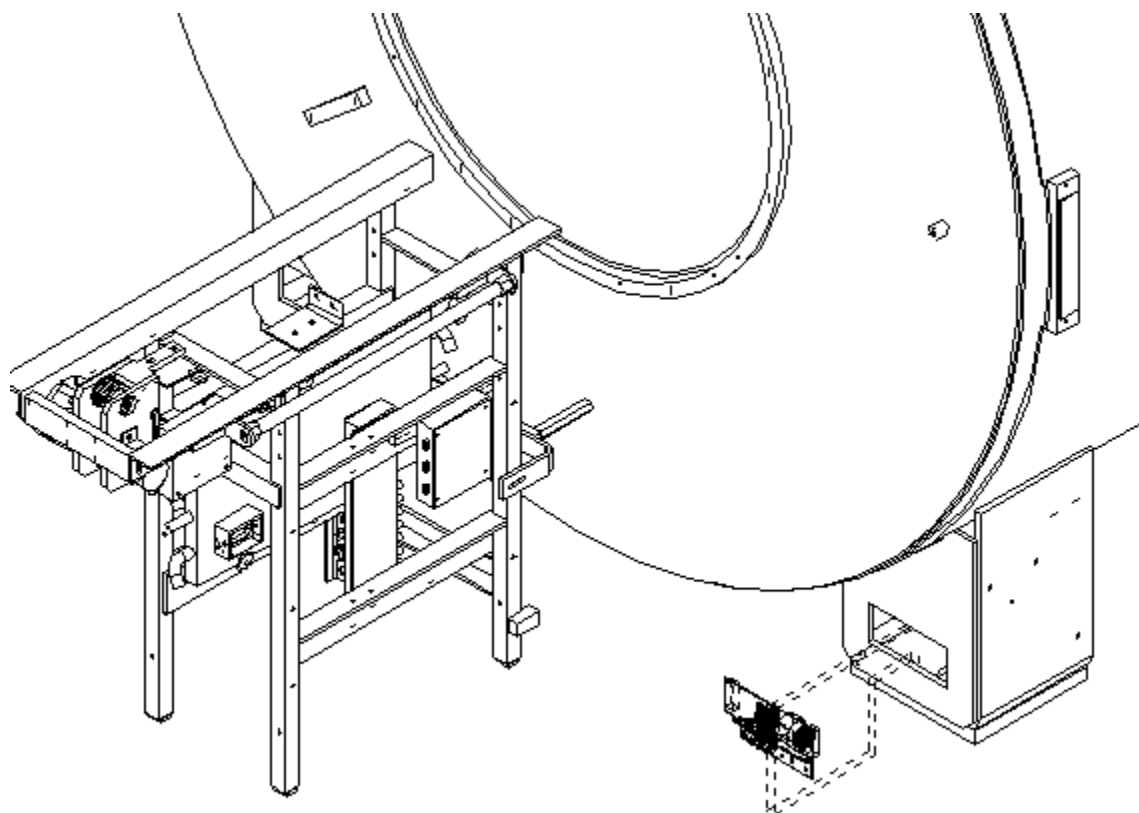

LIGHT SOURCE
ILLUSTRATION 1-2


8. Remove the comfort module by removing the two nuts at the top and the clips on the end bell. See Illustration 1-3. Disconnect the speaker cable that's inside of the comfort module. The light source cable will stay with the lights, which are attached to the RF coil.

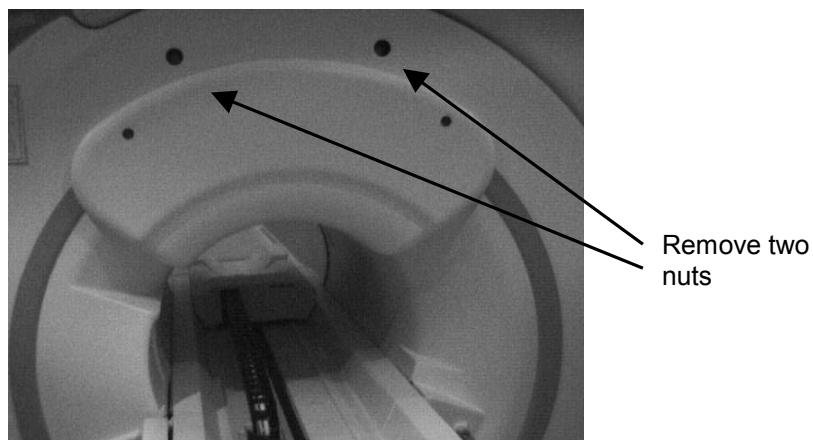

COMFORT MODULE
ILLUSTRATION 1-3


9. Slide rear pedestal back 2 feet.

10. Extend the RF Tube feet down all the way at the rear end bell by loosening the nut, and then extending the screws. The screws should be extended until the foot touches the TRM coil. See Illustration 1-4.

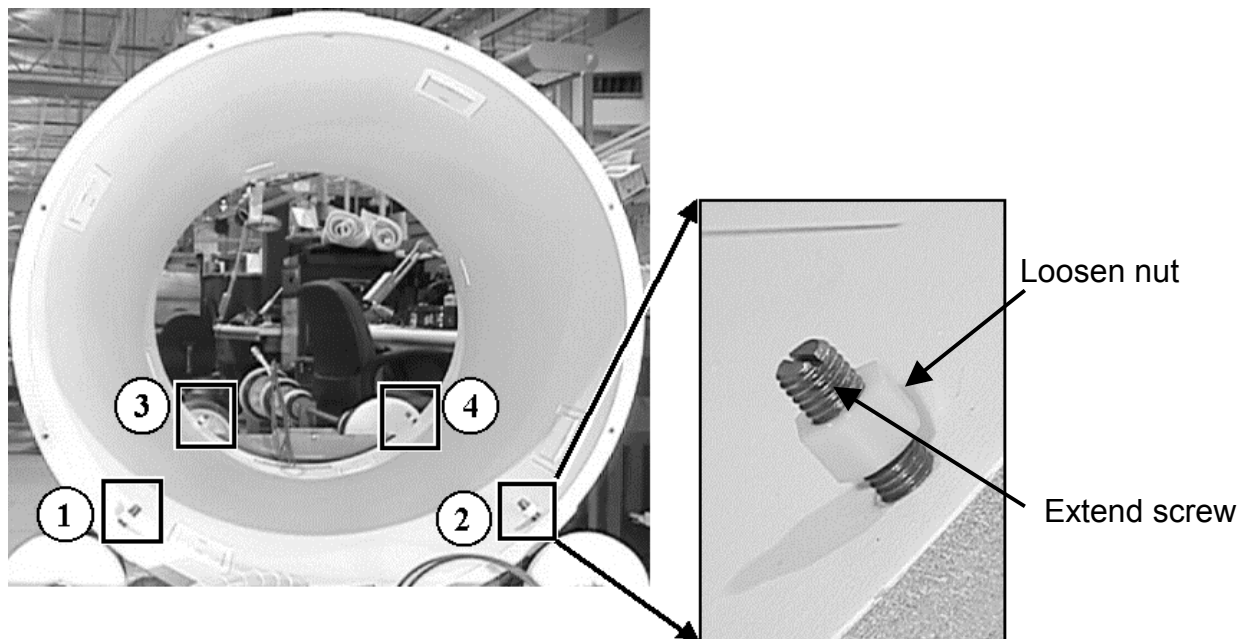

COIL FEET EXTENSION
ILLUSTRATION 1-5

11. Disconnect the vacuum hose at the rear end bell. See Illustration 1-6.

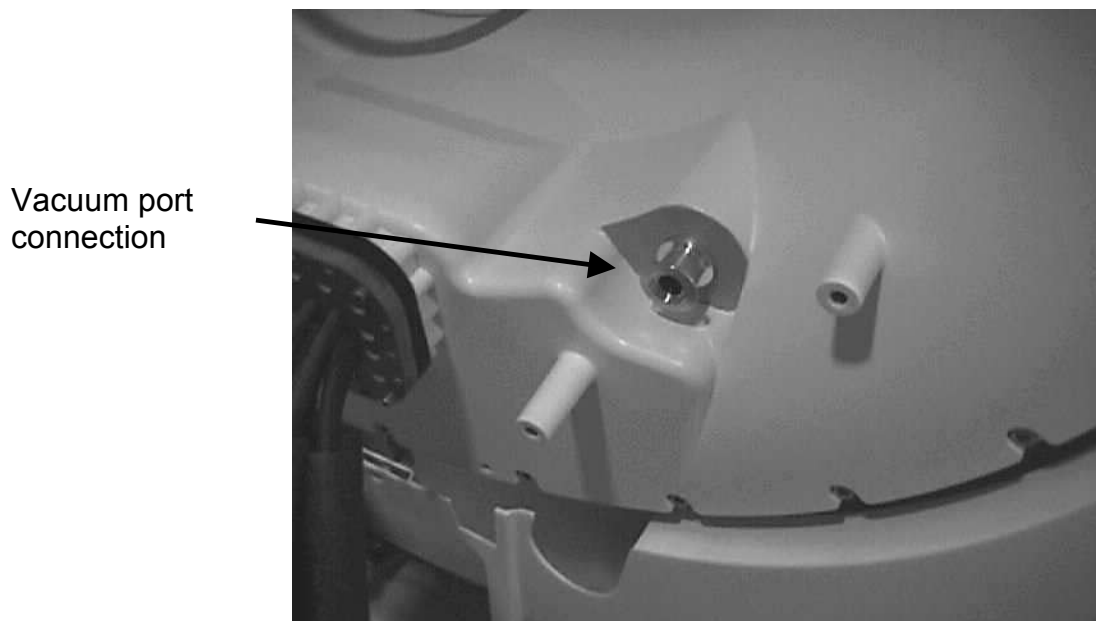


REAR END BELL
ILLUSTRATION 1-6

12. If you are replacing the end bell or removing the gradient coil, cut the 12 Gradient cables at the magnet room terminal block. Disconnect the optional resistive shim cable if

necessary.

13. Remove the 32 screws from the Gradient Vacuum seal fixture. See Illustration 1-7. Slide the cover plate back.

Location of 32
screws for gradient
vacuum seal

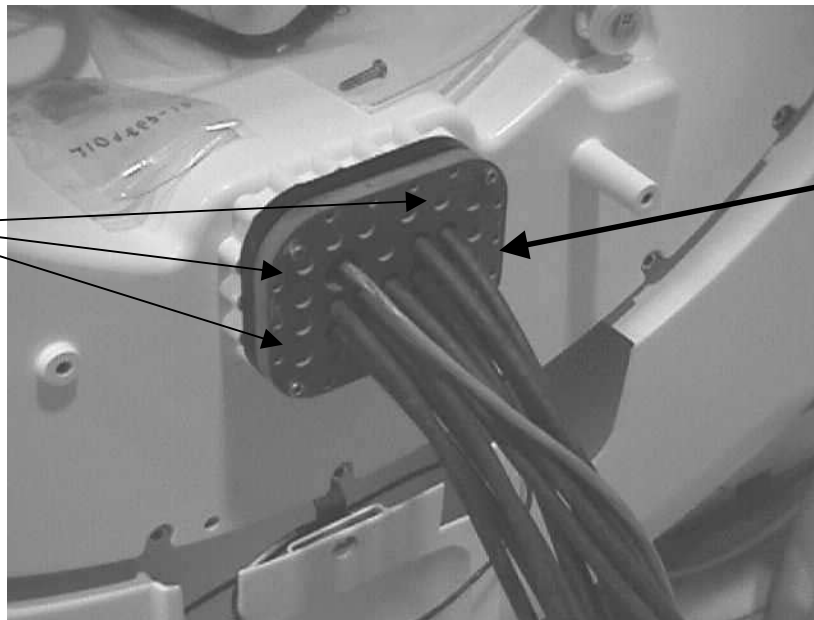

Cover plate

GRADIENT VACUUM SEAL FIXTURE
ILLUSTRATION 1-7

14. Remove the 8 nuts that secure the rear end bell to the RF Coil. See Illustration 1-8.

End bell



RF coil

NUTS SECURING
RF COIL TO END
BELL

REAR END BELL
ILLUSTRATION 1-8

15. Use an 8 mm allen wrench to remove the 24 M10 bolts that secure the rear end bell to the magnet.
16. Carefully pull the rear end bell back.

2- REAR END BELL INSTALLATION

1. Route the 12 Gradient cables and resistive shim cable through the opening in the rear end bell. Install all gaskets that were removed prior to terminating the gradient cables
2. Slowly move the rear end bell into place.
3. Use an 8 mm allen wrench to install the 24 M10 bolts that secure the rear end bell to the magnet. Use a star pattern when tightening the bolts. The front enclosure should touch the aluminum ring when finished.
4. Install the 8 nuts that secure the rear end bell to the RF Coil. See Illustration 1-8.
5. Install the 32 screws from the Gradient Vacuum seal fixture. See Illustration 1-7..
6. If necessary connect the optional resistive shim cable.
7. Install the vacuum pump connector at the rear end bell and install the cosmetic cover. See Illustration 1-6.
8. Re-terminate the gradient cables. Refer to Direction 2286976 for installation of Compression Splice and Cold Shrink.
9. Slide the rear pedestal to its proper location in back of the magnet.
10. Install the Patient Comfort Module.
11. Route and connect the cable going from the patient lights to the light source. See Illustration 1-2.
12. Connect the speaker and microphone cables to the intercom board located on rear pedestal. See Illustration 1-1.
13. Install the bridge. See procedure for Bridge Removal Twin Enclosure (ME2REA13.DOC).
14. Remove the lock out and tag out.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	Oct. 22, 2001	K. Schraufnagel	Preliminary Release
1	Nov 20, 2001	P. Kargard	Updates from Engineering