

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- OVERVIEW.....	2
2- BRIDGE REMOVAL	2
REVISION HISTORY	4

1- OVERVIEW

This procedure describes the replacement of the Temperature Sensor located in the front end bell of the open enclosure for Signa systems with the LCC Wide Open Magnet.

FRU	Part Number
Temperature Sensor Cable Kit	46-328326G1
Kit Contents:	
Run 749, MG2-A33-J10 TO MG2-A36, 35'	46-328354P1
Adhesive-Backed Fastener, 2" x 3"	46-328325P1
Self-Locking, Cable Tie 7.31" x .184"	46-208758P3

Tools Required:


Digital Volt Meter (DVM)
Standard Screwdriver
Philips Screwdriver
Adjustable Hex Wrench or Socket Set
Allen Wrench Set
2" x 4" piece of wood (approximately 18 - 20 inches long)

Manpower Requirements:

Two to three people required for removing end bells.

2- BRIDGE REMOVAL

1. Undock and move Patient Transport out of the way of the front cover of the magnet.
2. While holding in the patient transport dock plunger (see Illustration 2-1), move the carriage to the rear of the magnet by pressing the **[In Fast]** button on the enclosure.

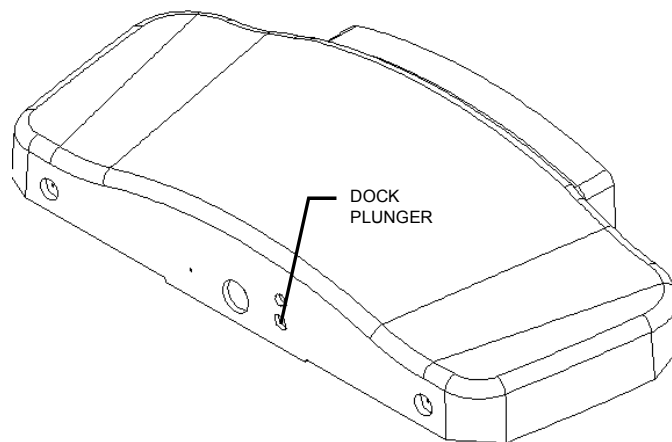

HOLDING DOCK PLUNGER
ILLUSTRATION 2-1

3. Lock out and tag out circuit breaker MR1 A14 CB6 at rear of RF/Pen Cabinet (MR1) on magnet enclosure power supply (MEPS) module.

For **RF/Pen II** and **RF/PDU Cabinet**: Lock out and tag out circuit breaker MR1 A20 CB4 at rear of RF/Pen II and RF/PDU Cabinet on System Support Module (SSM).

4. Using the DVM, verify that voltage between MR1 A14 TP1 (+) and TP14 (–) is 0 Vdc.

For **RF/Pen II** and **RF/PDU Cabinet**, Verify voltage between MR1 A20 TP1 (+) (Dock) and TP14 (–) (Dock RTN) is 0 Vdc.


5. Remove the bridge. See Bridge Replacement Twin Enclosure, ME2REA13.DOC.

6. Remove the cosmetic cover between the front end bell and the RF Coil by removing one of the Patient Light supports. The temperature sensor is located on the upper right side of the RF Coil.

7. Remove the tape that secures the temperature sensor to the RF Coil.

8. Remove the temperature sensor and the cable from the front end bell.

9. Install the replacement sensors into the holes in the RF Coil. See Illustration 2-2.

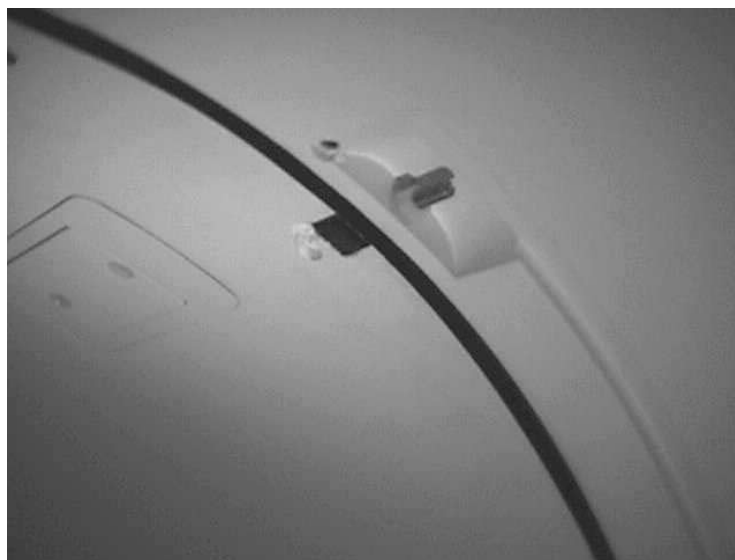

RF COIL
ILLUSTRATION 2-2

10. Install the provided tape over the temperature sensor.

11. Route the cable into the groove on the front end bell and reconnect the cable.

12. Re-install the Patient Light support, then re-install the cosmetic cover.

13. Re-install the bridge. See Bridge Replacement Twin Enclosure, ME2REA13.DOC.

14. Remove lock out and tag out circuit breaker MR1 A14 CB6 at rear of RF/Pen Cabinet (MR1) on magnet enclosure power supply (MEPS) module.

For **RF/Pen II** and **RF/PDU Cabinet**: Lock out and tag out circuit breaker MR1 A20 CB4 at rear of RF/Pen II and RF/PDU Cabinet on System Support Module (SSM).

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	Oct. 22, 2001	K. Schraufnagel	Preliminary Release