

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- LASER LIGHT CHECKS	2
REVISION HISTORY	4

1- LASER LIGHT CHECKS

Tools Required


- SPT Short Loader Cylinder, 2135652

Note


SPT Cylinder Availability - The SPT Short Loader Cylinder is part of the SPT Short Loader Assembly (2125244). This hardware, though shipped with Signa systems, is the property of GE. It is available for customer purchase.

The three laser lights, mounted on the front face of the magnet bore, are preadjusted at the factory and do not normally require realignment after the initial installation alignment procedures are completed.

1. Place SPT Short Loader Cylinder on patient table.
2. Press Align On button on magnet enclosure.
3. In the event that the laser beams are not aligned as shown in Illustration 1-1 and Illustration 1-2, refer to procedure for Patient Alignment Lights, Section 1- Laser Light Alignment Procedures.

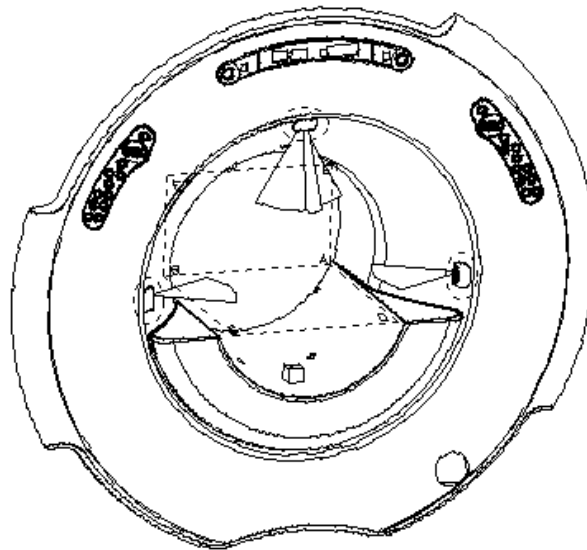

CHECKING ALIGNMENT LIGHTS
ILLUSTRATION 1-1

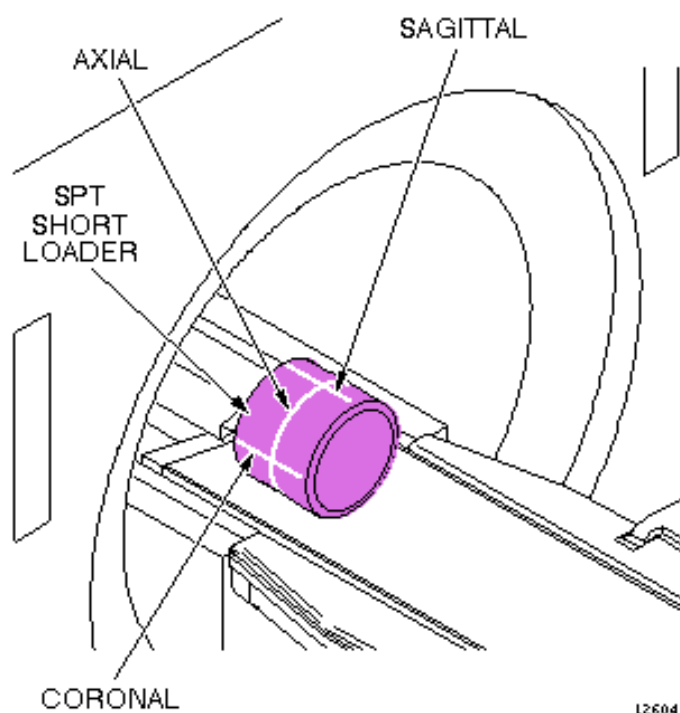

WRAP-AROUND BEAM PATTERN
ILLUSTRATION 1-2

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
A	Nov 1, 1998	W. Pasciak	Initial revision.
0	Feb 4, 1999	W. Pasciak	Final release