GE Healthcare

Product Network Filters (PNF)

IRIX version 3.0

GE CONFIDENTIAL

© 2006 General Electric Company

CLASS M DOCUMENTATION
PROPERTY OF GENERAL ELECTRIC COMPANY
FOR USE BY GE PERSONNEL ONLY
ANY OTHER USE EXPRESSLY PROHIBITED

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LANGUAGE

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

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This page left intentionally blank.

Revision History

Rev	Date	Reason For Change
1	12/07/2006	Initial Release

This page left intentionally blank.

Table of Contents

1 Product Network Filters for IRIX 11

 1.1 Overview 11

 1.2 Firewall State 12

 1.3 Firewall Settings 12

 1.4 Tabs 14

 1.4.1 Named Services Tab 14

 1.4.2 Allowed Nodes Tab 16

 1.4.3 DICOM Tab 18

 1.4.4 Expert Tab 19

 1.5 GE Healthcare Service Key Enabled View 21

 1.6 Troubleshooting 22

 1.7 Finalization 23

This page left intentionally blank.

1 Product Network Filters for IRIX


Last Revised:	12/07/2006
---------------	------------

1.1 Overview

Product Network Filters (PNF) is a software firewall tool which allows a site Network Administrator to control traffic from the site's internal network to the Diagnostic Imaging System. This tool was created in response to the customer's need for security.

This tool is accessible from the Common Service Desktop (CSD.)

Illustration 1: Main PNF Screen

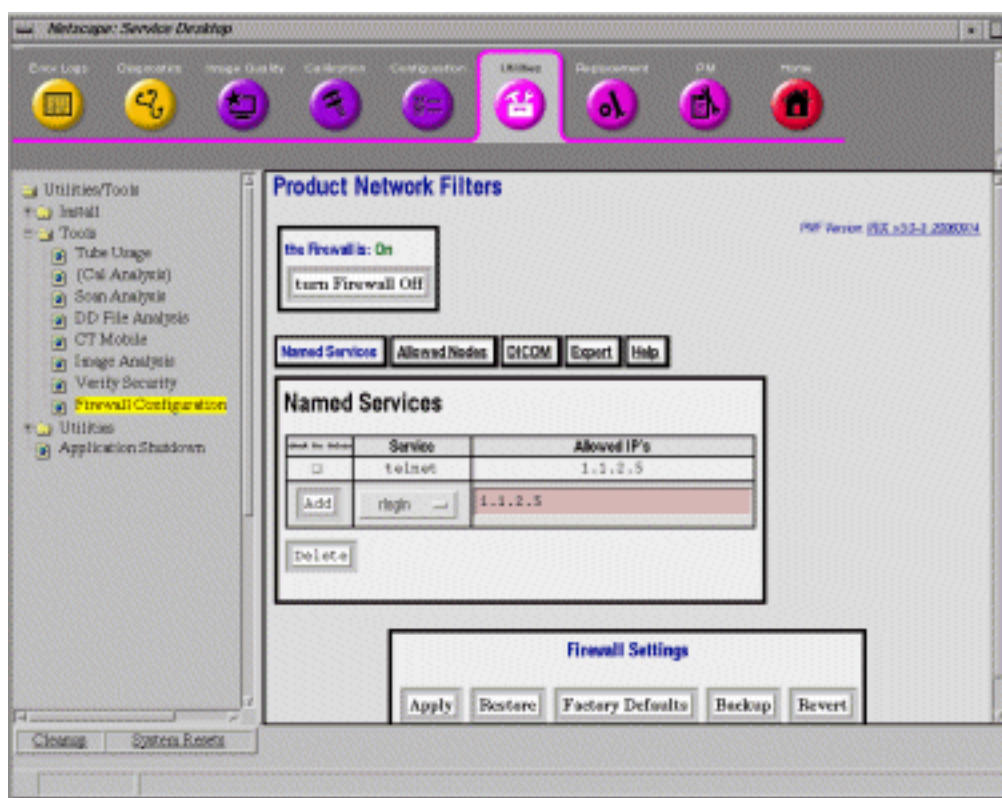The PNF offers a feature for GE use ONLY with the Insite connectivity configuration. The default value only allows the connection of the GE Healthcare Online Center (OLC) back office.

NOTE: Some customers have custom infrastructures and will want to change the IP address from where the service communication comes from. In this situation, we need to know from the customer these new addresses and configure them in the *GE Service* tab of PNF for remote service to work properly. The OLC connectivity team should be consulted in these situations.

**NOTICE**

ALL settings made in the PNF Tool are saved in the system save data. If a load from cold should need to be performed, the PNF settings will be restored during the system restore.

1.2 Firewall State

**NOTICE**

The Firewall should be ON at all times, this is the default setting. If issues are found because of the Firewall, contact the OLC connectivity team to resolve issues. It is not acceptable to leave the firewall OFF.

The Firewall can be turned ON and OFF by clicking the toggle button in [Illustration 2](#).

Illustration 2: Firewall State

1.3 Firewall Settings

The *Firewall Settings* toggle buttons are common across all tabs in the PNF tool. See [Illustration 3](#). These buttons can be clicked at any time and in any tab. The function of that button will be applied by the tool at the current state of **ALL** tabs.

Illustration 3: Firewall Settings Buttons

If there have not been any changes made the [Apply] and [Revert] buttons will not be active. The buttons will then appear as shown in [Illustration 4](#).

Illustration 4: Firewall Settings Buttons Inactive**Apply**

The PNF Tool saves the current filter settings (from all Tabs). If the firewall is currently "On", the application will restart the firewall with these new settings. The [Apply] button appears flat & grayed out when PNF detects that the current GUI settings are the same as those last applied.

Restore

The PNF Tool changes the GUI filter settings to those last saved with the [Backup] button. If the firewall is currently "On", the application will restart the firewall with these new settings.

Factory Defaults

The [Factory Defaults] button will reset the PNF configuration defaults. If the firewall is currently "On", the application will restart the firewall with these new settings.

Backup

The PNF Tool saves a copy of the filtering rules that were last applied via the [Apply] button. The firewall is not restarted.

Revert

The PNF Tool changes the GUI filter settings to those last made when the [Apply] button was clicked. The firewall is not restarted. The [Revert] button appears flat & grayed out when PNF detects that the current GUI settings are the same as those last applied.


1.4 Tabs

There are two modes (customer and GE service). In the service mode, the user can see all of the customer tabs as well.

1.4.1 Named Services Tab

The *Named Services* Tab, is the screen where the site's network administrator selects a standard networking service to open to this system from the site's network. See [Illustration 5](#).

Illustration 5: Named Services Tab

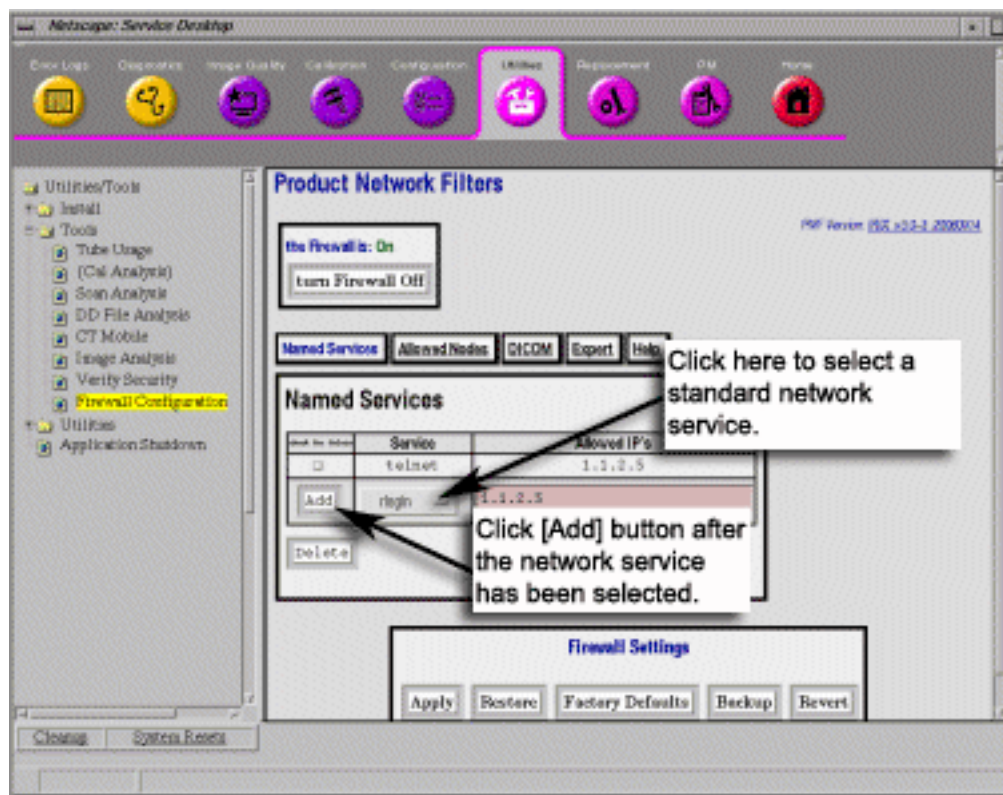Adding a Network Service

1. To add a network service, click the drop-down menu labeled *Services*. Standard Networking Services offered are:
 - telnet
 - ssh
 - ftp
 - rexec
 - rlogin
 - rsh
 - http
 - https
 - tftp

- ldap
 - ssl/tls
 - ntp
 - portmap
 - ping
2. Select a network service from the drop-down menu.
 3. An IP address can be entered in the *Allowed IP's* field or the field can be left blank. IP addresses allow the user to create rules where network communication is only accepted from a predetermined host or set of hosts. Below is a list of acceptable ways to populate the IP fields. The default value for IP fields is all IP addresses.
 - w.x.y.z = single IP address (ex: 1.2.34.5)
 - w.x.y.z/a.b.c.d = subnet with netmask (ex: 1.2.0.0/255.255.0.0)
 - w.x.y.z/a = subnet with byte netmask, this is the same as the previous subnet but specifying the netmask as a single byte (ex: 1.2.0.0/16)
 - Unlike PNF-Linux, the following syntax is not supported:
w.x.y.a - w.x.y.b = IP range, all hosts between two end points (ex: 1.2.34.5 - 1.2.34.10)
 4. Click [Add]. This service is now added to the list. If the *Allowed IP's* field was left blank an * will appear.
 5. Once all rules have been made, click the [Apply] to make changes to the PNF configuration.


Removing a Network Service

1. Check the box next to the *Service* on the list.
2. Click the [Delete] button.
3. Once all rules have been made, click the [Apply] to make changes to the PNF configuration.

1.4.2 Allowed Nodes Tab

The *Allowed Nodes* tab opens all network services to a specific IP address. See [Illustration 6](#).

Illustration 6: Allowed Nodes

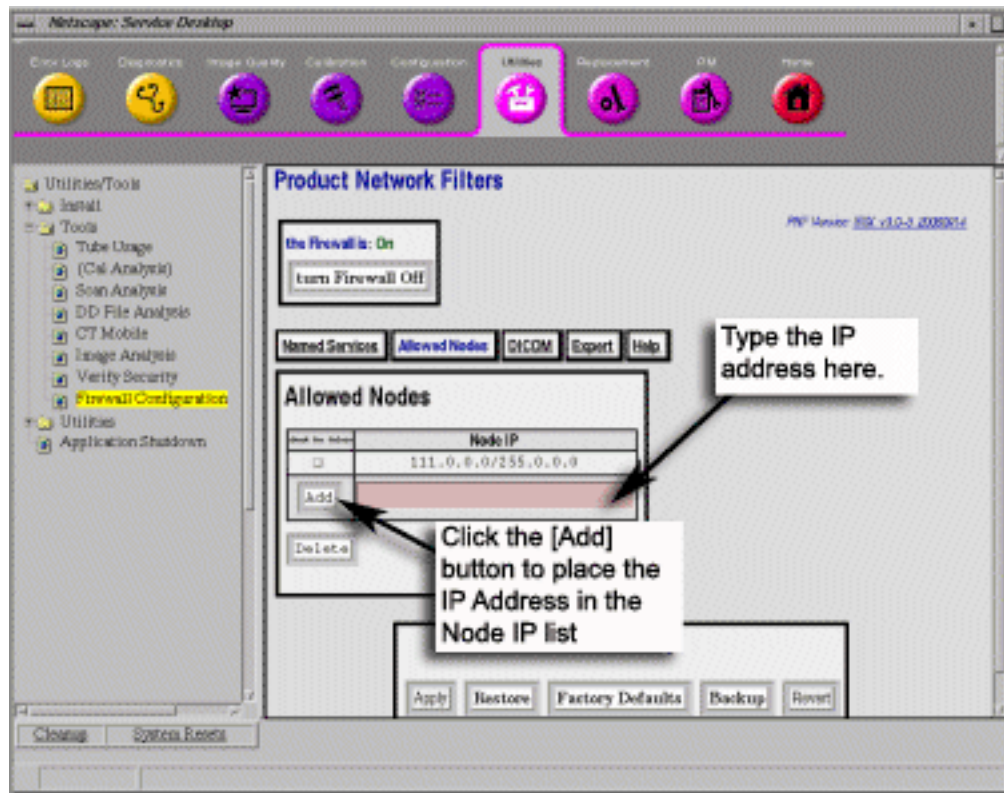Adding an Allowed Node

1. Enter the IP address in the *Node IP* field.
IP addresses allow the user to create rules where network communication is only accepted from a predetermined host or set of hosts. Below is a list of acceptable ways to populate the IP fields. The default value for IP fields is all IP addresses.
 - w.x.y.z = single IP address (ex: 1.2.34.5)
 - w.x.y.z/a.b.c.d = subnet with netmask (ex: 1.2.0.0/255.255.0.0)
 - w.x.y.z/a = subnet with byte netmask, this is the same as the previous subnet but specifying the netmask as a single byte (ex: 1.2.0.0/16)
2. Once the IP address has been entered, click [Add] so the IP Address will be added to the list.
3. Once all Nodes have been added, click the [Apply] button to make changes to the PNF configuration.

Removing an Allowed Node


1. Check the box next to the *Node IP* in the list.
2. Click the [Delete] button to remove the *Node IP* from the list.

3. Once all changes have been made, click the [Apply] button to make changes to the PNF configuration.

1.4.3 DICOM Tab

This screen may be used to verify that the correct DICOM port was populated during installation. This field should be auto-populated during the initial set-up of the DICOM ports. Also additional DICOM ports can be added. See [Illustration 7](#).

Illustration 7: DICOM Tab

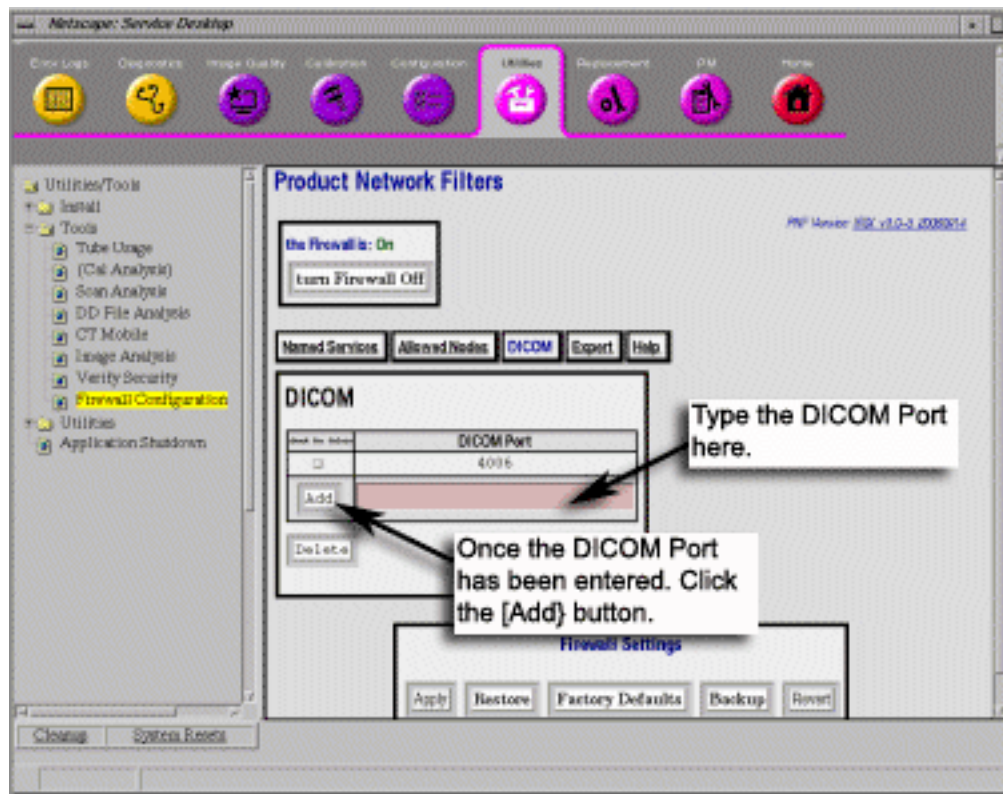Adding a DICOM Port

1. Enter the DICOM Port in the field next to the [Add] button.
2. Click the [Add] button so the new DICOM Port is in the list.
3. Once all DICOM Ports have been added, click the [Apply] button to make changes to the PNF configuration.


Removing a DICOM Port

1. Check the box next to the *DICOM Port* in the list.
2. Click the [Delete] button to remove the *DICOM Port* from the list.
3. Once all changes have been made, click the [Apply] to make changes to the PNF configuration.

1.4.4 Expert Tab

The *Expert* Tab allows a more flexible method for the network administrator to create rules. Any input to this tab will require networking expertise, please work with the site's network administrator when working in this tab.

Illustration 8: Expert Tab

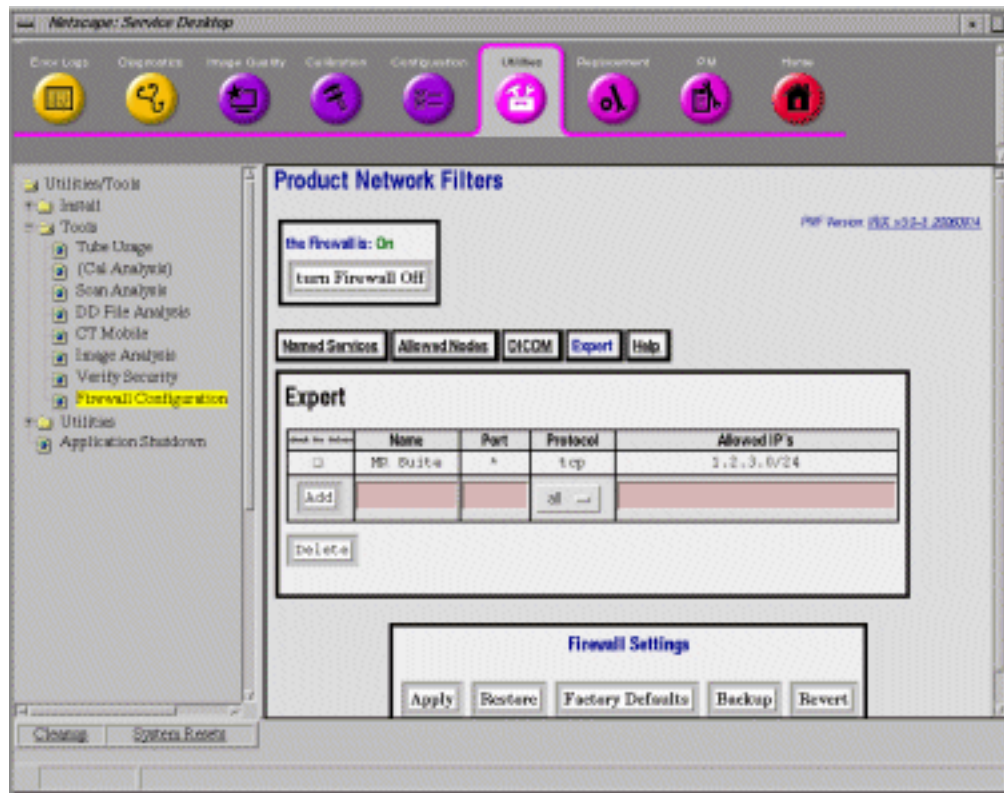Creating a Rule

1. To create a rule for the firewall, enter the application name in the *Name* field.
2. Enter the Port info in the *Port* field.
3. Click on the *Protocol* drop-down menu for the rule to be applied.
4. If needed, enter the IP address that this rule will be applied.
IP addresses allow the user to create rules where network communication is only accepted from a predetermined host or set of hosts. Below is a list of acceptable ways to populate the IP fields. The default value for IP fields is all IP addresses.
 - w.x.y.z = single IP address (ex: 1.2.34.5)
 - w.x.y.z/a.b.c.d = subnet with netmask (ex: 1.2.0.0/255.255.0.0)
 - w.x.y.z/a = subnet with byte netmask, this is the same as the previous subnet but specifying the netmask as a single byte (ex: 1.2.0.0/16)
5. Once all fields have been populated, click the [Add] button. The rules created will now be added to the list.

6. Once all rules have been made, click the [Apply] to make changes to the PNF configuration.

Removing a Rule


1. Check the box next to the *Name* in the list.
2. Click the [Delete] button to remove the rule from the list.
3. Once all changes have been made, click the [Apply] to make changes to the PNF configuration.

1.5 GE Healthcare Service Key Enabled View

The *GE Service* Tab displays the IP address for the connection of remote service. The ONLY time this address would be changed is if the site is using a non-standard VPN administrator. IF this is the case, obtain the correct IP address from the site network administrator. Contact the OLC Connectivity team to properly establish the link. See [Illustration 9](#).

The PNF offers a feature for GE use ONLY with the Insite connectivity configuration. The default value only allows the connection of the GE Healthcare Online Center (OLC) back office.

Illustration 9: GE Service Tab

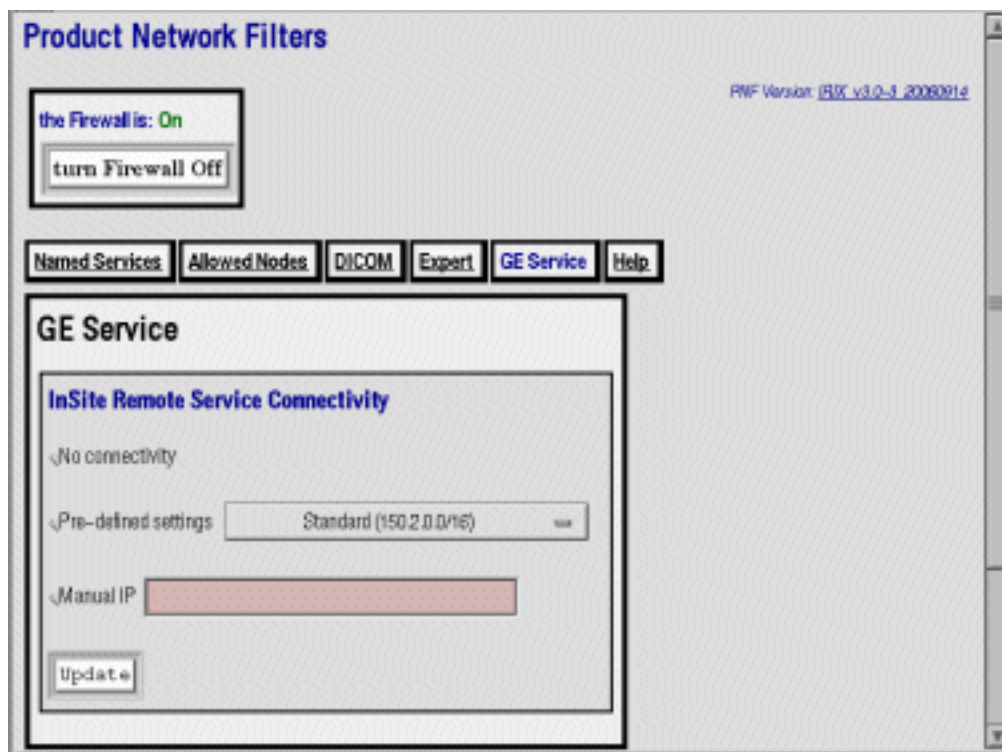1.6 Troubleshooting

If the Firewall is suspected of causing an issue. A method to check if the Firewall is the source, turn the Firewall OFF and see if the issue remains. Make sure that the Firewall is ON before exiting. See [Illustration 10](#) .

Illustration 10: Firewall State Button

The Firewall can be turned ON and OFF by clicking the [turn Firewall Off] or [turn Firewall On] buttons.

NOTICE

The Firewall should be ON at all times, this is the default setting. If issues are found because of the Firewall, contact the OLC connectivity team to resolve issues. It is not acceptable to leave the firewall OFF.

1.7 Finalization

Once all the necessary changes/settings have been made, perform the following steps:

1. At the *Firewall Setting* buttons, click [Apply]. See [Illustration 11](#) . This will apply the changes to the PNF configuration.
2. In order to store this configuration in the system data, click the [Backup] button. See [Illustration 11](#) . If a load from cold should need to be performed, the PNF configuration settings will be restored from the SaveInfo media.

Illustration 11: Firewall Setting Buttons

This page left intentionally blank.

© 2006, General Electric Company.
GE Medical Systems, a General Electric Company, going to market as GE Healthcare.
3200 N. Grandview Blvd
Waukesha, WI 53186
US

www.gehealthcare.com