

TABLE OF CONTENTS


TABLE OF CONTENTS	1
1- INTRODUCTION	2
2- OPERATOR WORKSPACE TABLE	3
2-1 Monitor Replacement	3
2-1-1 Tools Required	4
2-1-2 Removal Procedures	4
2-1-3 Installation Procedures	4
2-2 Workspace Interface Module Replacement	5
2-2-1 Tools Required	5
2-2-2 Removal Procedures	5
2-2-3 Installation Procedures	6
2-3 DAT Drive Replacement	7
2-3-1 Tools Required	7
2-3-2 Removal Procedures	7
2-3-2 Removal Procedures	8
2-3-3 Installation Procedures	8
2-4 DAT Legacy Enclosure Replacement	9
2-4-1 Tools Required	9
2-4-2 Removal Procedures	9
2-4-3 Installation Procedures	10
2-5 Pioneer MOD Drive Replacement	11
2-5-1 Tools Required	11
2-5-2 Removal Procedures	11
2-5-3 Installation Procedures	12
2-6 Pioneer MOD Legacy Enclosure Replacement	13
2-6-1 Tools Required	13
2-6-2 Removal Procedures	13
2-6-3 Installation Procedures	14
2-7 Modem Replacement	15
2-7-1 Tools Required	15
2-7-2 Removal Procedures	15
2-7-3 Installation Procedures	16
2-8 Digital DASM Replacement	17
2-8-1 Tools Required	17
2-8-2 Removal Procedures	17
2-8-3 Installation Procedures	18
2-9 Analog DASM Replacement	19
2-9-1 Tools Required	19
2-9-2 Removal Procedures	19
2-9-3 Installation Procedures	20
2-10 LCD Panel Replacement	21
2-10-1 Tools Required	21
2-10-2 Removal Procedures	21
2-10-3 Installation Procedures	22
2-11 Keyboard / Mouse Replacement	23
2-11-1 Tools Required	23
2-11-2 Removal Procedures	23
2-11-3 Installation Procedures	24
REVISION HISTORY	25

DESCRIPTION

This procedure describes the replacement procedures for all FRU-1 parts for the Operator Workspace. For replacement of 8.X Indigo computer parts, refer to procedure "Replacement and Maintenance 8X Indigo Computer."


1- INTRODUCTION

The operator workspace (OW) (see Illustration 1-1) subsystem is serviced by isolating faults to a FRU level, and by replacing the defective FRU. This procedure describes the replacement of level 1 FRUs. Section 2 covers procedures for items associated with the OW Table, while section 3 covers procedures for items associated with the OW Cabinet.

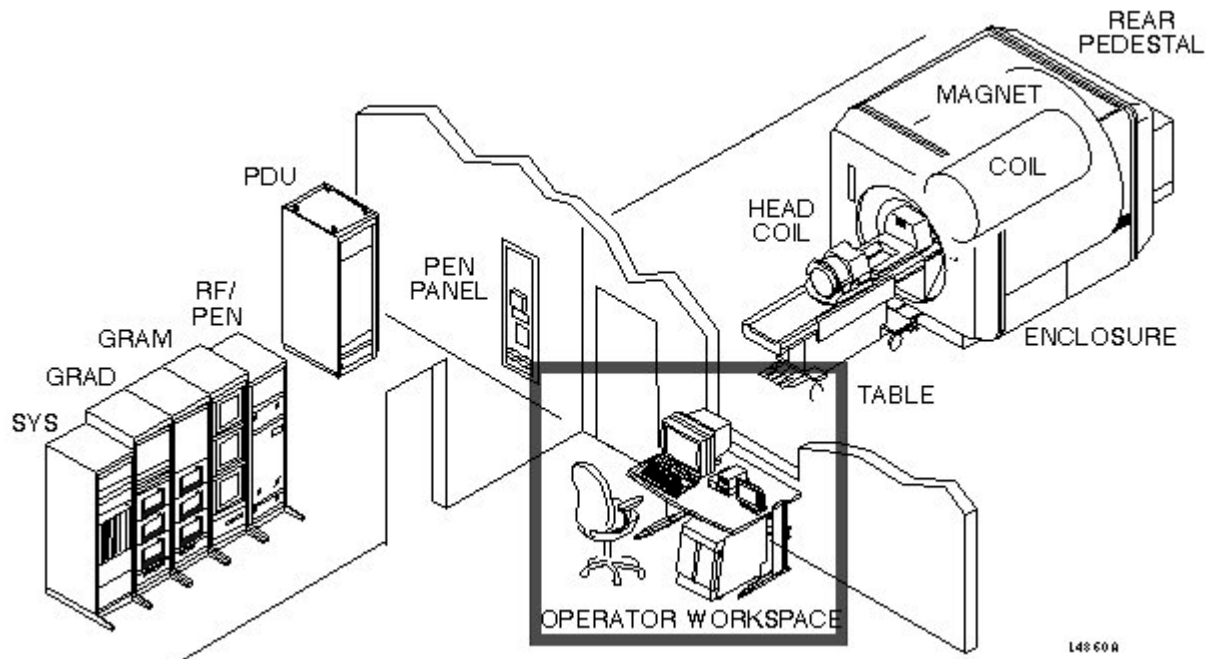

OPERATOR WORKSPACE
ILLUSTRATION 1-1

2- OPERATOR WORKSPACE TABLE

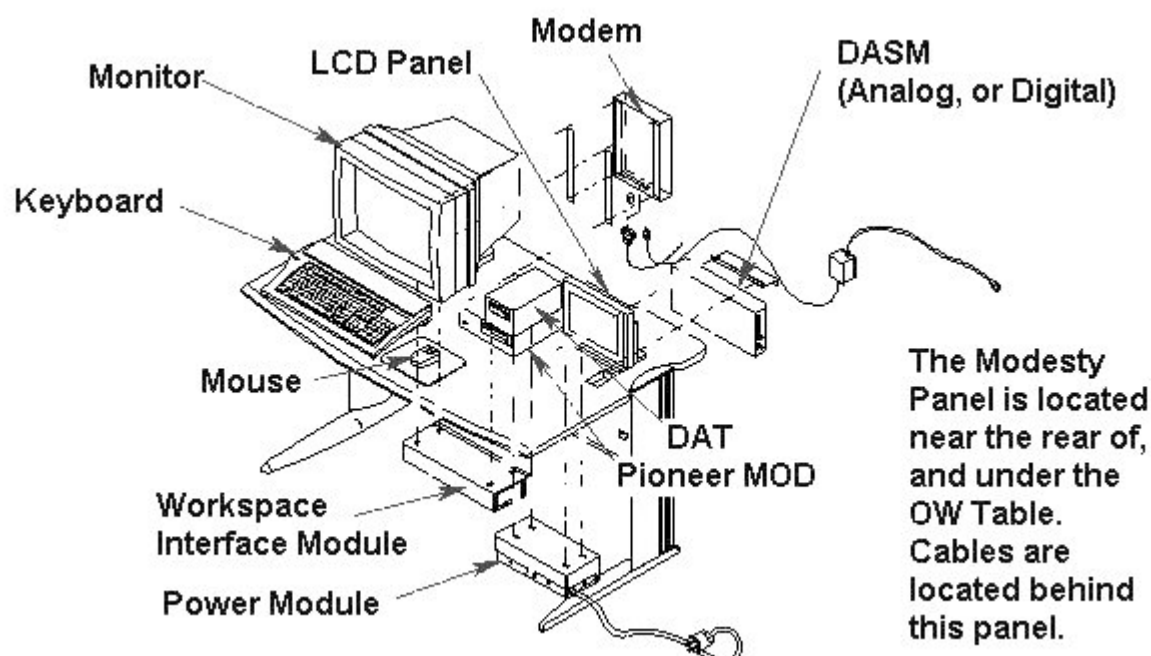

OPERATOR WORKSPACE TABLE
ILLUSTRATION 2-1

2-1 Monitor Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Log out from Host Computer.
2. Shut down the Host Computer.
3. Bring the PC down by closing the Program Manager. Be sure that no LEDs are lighted.
4. As applicable, open the circuit breaker(s) that power the Operator Workspace to turn power off.
5. Lock and tag the circuit breaker(s). (Refer to CD-ROM *Dir. 2187583-3 [or -2], MR Release Signa 5x/8x Service Methods, Renewal Parts and Service Tools*, Safety, Section 6, OSHA LOCKOUT/TAGOUT REQUIREMENTS.)
6. Open computer cabinet door and verify that power is off by twice pressing (on/off) the Host Computer power button; the LED should remain off.
7. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.

2-1-1 Tools Required

- Small flat head screwdriver
- Medium flat head screwdriver

2-1-2 Removal Procedures

1. Locate the monitor on the left side of the workspace table. See Illustration 2-1 (page 3 of this document).
2. At the rear of the monitor, disconnect the external cables.
3. Remove the old monitor from the workspace table and set aside.

2-1-3 Installation Procedures

1. Remove the replacement Monitor from the ESD Storage Bag and/or original carton, as applicable.
2. Securely position the Monitor on the left side of the workspace table.
3. Reconnect the external cables to the monitor.
4. See Procedure for Monitor Adjustments (Gray Scale) or Procedure for Monitor Adjustments (Color) for monitor setups.
5. Return power to the unit. Place the old monitor in the ESD Storage Bag and/or original carton, as applicable.

2-2 Workspace Interface Module Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-2-1 Tools Required

- Small flat head screwdriver
- Medium flat head screwdriver

2-2-2 Removal Procedures

1. Locate the Workspace Interface Module (WIM) on the left underside of the workspace table. See Illustration 2-1 (on page 3 of this document).
2. From the front of the WIM, disconnect the keyboard cable.
3. From the rear of the WIM, disconnect the following cables: J7 Patient Interface, J8 Host Interface, J9 PC Interface, J17 Power, and J18 TYME I/F. If present, remove all other cables residing in ports J13 through J16.

2-2-2 Removal Procedures (continued)

4. Remove the WIM by turning each screw on the rear of the Interface clockwise 1/2 turn. With the screw held in a horizontal position, pull the module toward you.



Possible equipment damage. Screws are spring loaded. Be sure to hold on to the module tightly as you turn the second screw. Failure to do so will cause the module to fall.

Repeat this procedure for the second screw.

5. Set the old WIM aside for later reattachment of the cover. See section 2-2-3 Installation Procedures, step 2.

2-2-3 Installation Procedures

1. Remove the replacement WIM from the ESD storage bag and/or original carton, as applicable.
2. On the replacement WIM, unscrew the screws located at each corner of the cover. Remove the module cover and reinstall it on the old module.
3. With the cover removed, set the jumpers on the replacement as follows: JP1, JP2, JP5, and JP7 should be in the *out* position. Jumpers JP3, JP4 (Autovoice), JP6, and JP8 should all be in the *in* position.
4. Reinstall the replacement WIM under the workspace table by sliding it into the existing cover from the rear. Make sure the module clears the top of existing cover and that the cable ports are accessible from the rear of the table. The module snaps into place and the screws automatically return to a vertical position.
5. Reconnect the external cables to the WIM. Make sure that all cables match the corresponding ports. See step 3 of section 2-2-2 Removal Procedures.
6. Return power to the unit.
7. Store the old WIM in the ESD storage bag and/or original carton, as applicable.
8. See Procedure for SGI Diagnostics for the keyboard to test the installation.

2-3 DAT Drive Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-3-1 Tools Required

- Medium Phillips screwdriver
- Small flat head screwdriver
- Grounding wrist strap
- ESD Kit
- Container for screws

2-3-2 Removal Procedures


1. Locate the DAT Drive on the workspace table. See Illustration 2-1 (on page 3 of this document).
2. At the rear of the DAT Drive, disconnect the external cables.
3. Unscrew the screws located at the bottom corners of the drive.
4. Remove the DAT Drive Enclosure (cover) by sliding it off the rest of the case.
5. Unscrew only the middle screws located on either side of the drive bracket.
6. Lift the DAT Drive up and to the right to display the drive cables.

2-3-2 Removal Procedures

7. Carefully disconnect the wide ribbon cable by moving it gently from side to side until it is removed. Then remove the white tip cable (white, red and black wire) and the small black tip cable (blue, red, yellow and brown wire).
8. Place the old DAT Drive on the red ESD mat.

2-3-3 Installation Procedures

1. Remove the replacement DAT Drive from the ESD Storage Bag and/or original carton, as applicable and set Jumper and switches.

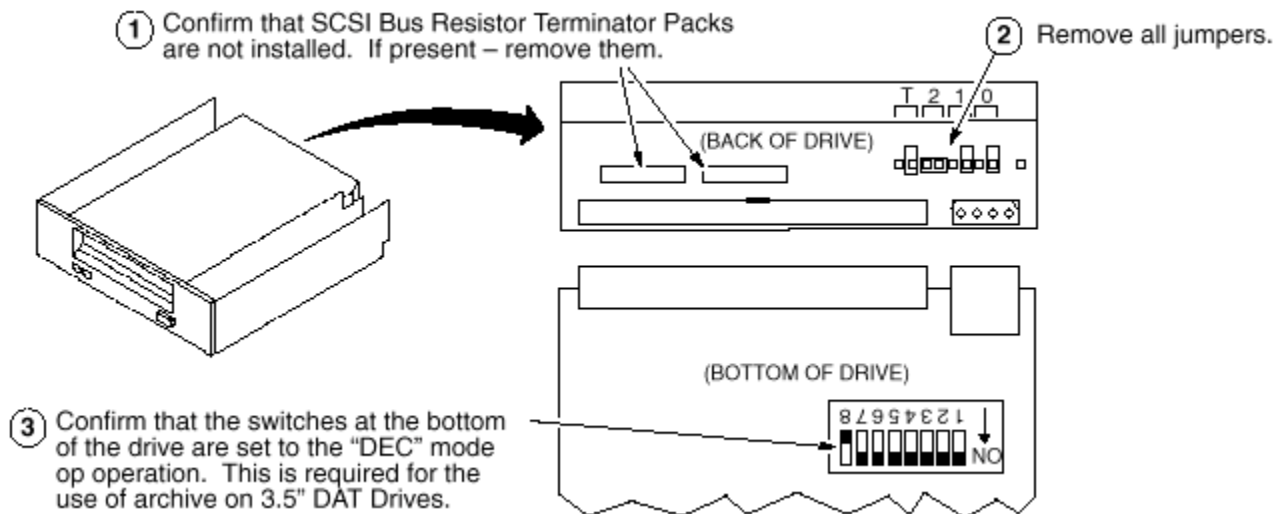

2. Reconnect the internal cables to the DAT Drive. Make sure the drive is facing up when you reconnect the cables. See step 7 of section 2-3-2 Removal Procedures.
3. Angle the DAT Drive down and to the left to place it back inside the case.
4. Reinstall the middle screws located on either side of the drive bracket.
5. Reinstall the DAT Drive Enclosure (cover) by sliding it over the rest of the case. Make sure the case vents appear at the top front of the drive to ensure proper alignment of the bottom screws.
6. Reinstall the bottom screws.
7. Position the DAT Drive in the original location.
8. Reconnect the external cables to the DAT Drive.
9. Return power to the unit.
10. Place the old DAT Drive in the ESD Storage Bag and/or original carton, as applicable.

2-4 DAT Legacy Enclosure Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-4-1 Tools Required

- Medium Phillips screwdriver
- Small flat head screwdriver
- Grounding wrist strap
- ESD Kit
- Container for screws

2-4-2 Removal Procedures

1. Locate the DAT Drive on the workspace table. See Illustration 2-1 (on page 3 of this document).
2. At the rear of the DAT Drive, disconnect the external cables.
3. Unscrew the screws on the bottom of the case to release the Drive Enclosure (cover).
4. Remove the old DAT Legacy Enclosure (cover) by sliding it off the rest of the case.
5. Place the old DAT Legacy Enclosure on the red ESD mat.

2-4-3 Installation Procedures

1. Remove the replacement DAT Legacy Enclosure from the ESD Storage Bag and/or original carton, as applicable.
2. Reinstall the DAT Legacy Enclosure (cover) by sliding it over the rest of the case. Make sure the case vents appear at the top front of the drive to ensure proper alignment of the bottom screws.
3. Reinstall the bottom screws.
4. Position the DAT Drive at the middle of the workspace table.
5. Reconnect the external cables to the DAT Drive.
6. Return power to the unit.
7. Place the old DAT Legacy Enclosure in the ESD Storage and/or original carton, as applicable.

2-5 Pioneer MOD Drive Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-5-1 Tools Required

- Small flat head screwdriver
- Medium Phillips screwdriver
- Grounding wrist strap
- ESD Kit
- Container for screws

2-5-2 Removal Procedures

1. Locate the Pioneer Drive on the workspace table. See Illustration 2-1 (on page 3 of this document).
2. At the rear of the Drive, disconnect the external cables.
3. Unscrew the screws on the bottom of the Drive that secure the case cover.
4. Remove the Drive Enclosure (cover) by sliding it away from the rest of the case.
5. To remove the Drive, first unscrew the bracket screws on the bottom of the case. Then lift the Drive up and out.
6. Unscrew the screws from the bracket to remove the Pioneer Drive.

2-5-2 Removal Procedures (continued)

7. Carefully disconnect the wide ribbon cable by moving it gently from side to side until it is removed. Then remove the white tip cable (white, red and black wire).

8. Place the old Pioneer Drive on the red ESD mat.

2-5-3 Installation Procedures

1. Remove the replacement Pioneer Drive from the ESD Storage Bag and/or original carton, as applicable.
2. Reconnect the internal cables to the Pioneer Drive. Make sure the drive is facing up when you reconnect the cables. See step 7 of section 2-5-2 Removal Procedures.
3. Position the Pioneer Drive on the bracket and reinstall the screws.
4. Reinstall the Pioneer Drive Enclosure (cover) by sliding it over the rest of the case. Make sure the case vents appear at the top front of the drive to ensure proper alignment of the bottom screws.
5. Reinstall the bottom screws.
6. Reconnect the external cables to the Pioneer Drive.
7. Position the Pioneer Drive at the middle of the workspace table.
8. Return power to the unit.
9. Place the old Pioneer Drive in the ESD Storage Bag and/or original carton, as applicable.

2-6 Pioneer MOD Legacy Enclosure Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-6-1 Tools Required

- Small flat head screwdriver
- Medium Phillips screwdriver
- Grounding wrist strap
- ESD Kit
- Container for screws

2-6-2 Removal Procedures

1. Locate the Pioneer Drive on the workspace table. See Illustration 2-1 (on page 3 of this document).
2. At the rear of the drive, disconnect the external cables.
3. Unscrew the screws on the bottom of the case to release the Pioneer Legacy Enclosure (cover).
4. Remove the Pioneer Legacy Enclosure (cover) by sliding it away from the rest of the case.
5. Place the old Pioneer Legacy Enclosure on the red ESD mat.

2-6-3 Installation Procedures

1. Remove the replacement Pioneer Legacy Enclosure ESD Storage Bag and/or original carton, as applicable.
2. Reinstall the Pioneer Legacy Enclosure (cover) by sliding it over the rest of the case. Make sure the case vents appear at the top front of the drive to ensure proper alignment of the bottom screws.
3. Reinstall the bottom screws.
4. Reconnect the external cables to the Pioneer Drive.
5. Position the Pioneer Drive at the middle of the workspace table.
6. Return power to the unit.
7. Place the old Pioneer Legacy Enclosure in the ESD Storage Bag and/or original carton, as applicable.
8. Reinstall system, power on and boot the system into **[root]**.
9. From the Root menu bar on the right hand side of the desktop, select **[Reconfigure]** by clicking the right mouse button over the reconfigure option.
10. Left mouse button click the option **[Configure Devices]**. A "Term" window will now open and remain blank for up to two minutes. When finished, "Done with Device Configuration" will appear in the Term window. Press any key to exit. Reboot the system to **[Signa]** and test the system.

2-7 Modem Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-7-1 Tools Required

- Small flat head screwdriver
- Medium Phillips screwdriver
- Container for screws

2-7-2 Removal Procedures

1. Locate the Modem at the right rear of the workspace table, attached to the modesty panel. See Illustration 2-1 (on page 3 of this document).
2. Release the Velcro straps to remove the Modem.
3. At the rear of the Modem, disconnect the external cables.
4. Remove the old Modem and set aside.

2-7-3 Installation Procedures

1. Remove the replacement Modem from the ESD Storage Bag and/or original carton, as applicable.
2. Reconnect the external cables to the Modem.
3. Reposition the Modem at the right rear of the workspace table, attached to the modesty panel.
4. Reattach the Velcro straps to secure the Modem.
5. Return power to the unit.
6. Place the old Modem in the ESD Storage Bag and/or original carton, as applicable.

2-8 Digital DASM Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-8-1 Tools Required

- Small flat head screwdriver
- Medium Phillips screwdriver
- ESD Kit
- Container for screws

2-8-2 Removal Procedures

1. Locate the Digital DASM on the right rear of the workspace table, in the cable tray. See Illustration 2-1 (on page 3 of this manual).
2. At the rear of the Digital DASM, disconnect the external cables.
3. Remove the Digital DASM.
4. Place the old Digital DASM on the red ESD mat.

2-8-3 Installation Procedures

1. Remove the replacement Digital DASM from the ESD Storage Bag and/or original carton, as applicable.
2. Reposition the Digital DASM on the right rear of the workspace table, in the cable tray.
3. Reconnect the external cables to the Digital DASM.
4. Return power to the unit.
5. Place the old Digital DASM in the ESD Storage Bag and/or original carton, as applicable.
6. Perform Digital DASM Calibration. See the procedure for the Digital DASM Calibration and Troubleshooting

2-9 Analog DASM Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6*.)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-9-1 Tools Required

- Small flat head screwdriver
- Medium Phillips screwdriver
- ESD Kit
- Container for screws

2-9-2 Removal Procedures

1. Locate the Analog DASM on the right rear of the workspace table, in the cable tray. See Illustration 2-1 (on page 3 of this document).
2. At the rear of the Analog DASM, disconnect the external cables.
3. Remove the Analog DASM.
4. Place the old Analog DASM on the red ESD mat.

2-9-3 Installation Procedures

1. Remove the replacement Analog DASM from the ESD Storage Bag and/or original carton, as applicable.
2. Reposition the Analog DASM on the right rear of the workspace table, in the cable tray.
3. Reconnect the external cables to the Analog DASM.
4. Return power to the unit.
5. Place the old Analog DASM in the ESD Storage Bag and/or original carton, as applicable.
6. Perform Analog DASM Calibration. See the procedure for the Analog DASM Calibration and Troubleshooting for setup information for this device.

2-10 LCD Panel Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6*.)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-10-1 Tools Required

- Medium Phillips screwdriver
- Small flat head screwdriver

2-10-2 Removal Procedures

1. Locate the LCD Panel on the right side of the workspace table. See Illustration 2-1 (on page 3 of this document).
2. At the rear of the LCD Panel, disconnect the external cables.
3. Remove the old LCD Panel and set aside.

2-10-3 Installation Procedures

1. Remove the replacement LCD Panel from the ESD Storage Bag and/or original carton, as applicable.
2. Reposition the replacement LCD Panel on the right side of the workspace table.
3. Reconnect the external cables to the LCD Panel.
4. Return power to the unit.
5. See Procedure for LCD Brightness/Contrast Adjustment to setup this device.
6. Place the old LCD Panel in the ESD Storage Bag and/or original carton, as applicable.

2-11 Keyboard / Mouse Replacement



FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE INVOLVED UNITS BEFORE STARTING THIS PROCEDURE.

1. Open computer cabinet door and verify that power to the Operator Workstation is off by twice pressing (on/off) the power button on the Host Computer; the power LED should remain unlighted. If not, use the lock out/tag out steps listed in Section 2-1, Monitor Replacement. (Refer to *Procedure For Safety: Section 6.*)
2. Use proper ESD precautions when working on this hardware. These components are ESD sensitive and damage will occur to the hardware if proper ESD measures are not followed.



Use caution when removing panels and modules from the Operator Workspace. ESD copper fingers located on the inside edges of the SGI housing and cover panels damage easily. Handle carefully to prevent damage to your hands and the ESD fingers. If damage occurs to the ESD fingers, they must be replaced in their slot like zipping a zipper. Start where the fingers are aligned correctly and press the damaged fingers in place by sliding your fingers along the edges as applicable.

2-11-1 Tools Required

- Small flat head screwdriver
- Medium Phillips screwdriver

2-11-2 Removal Procedures

1. Locate the Keyboard/Mouse on the workspace table. See Illustration 2-1 (on page 3 of this document).
2. Turn the Keyboard over to display the cable connections on the underside.
3. Disconnect the external cables.
4. Remove the old Keyboard and Mouse and set aside.

2-11-3 Installation Procedures

1. Remove the replacement Keyboard/Mouse from the carton and/or ESD Storage Bag as applicable.
2. Reposition the Keyboard/Mouse on the workspace table.
3. Reconnect the external cables to the Keyboard.
4. Return power to the unit.
5. See Procedure for SGI Diagnostics, keyboard diagnostics to test the installation.
6. Place the old Keyboard/Mouse in the ESD Storage Bag and/or original carton, as applicable.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
A	Mar 23, 1998	R. Hawthorne	Initial Conversion From ToolBook To Word
0	Apr 2, 1998	R. Hawthorne	Added system configuration to MOD Replacements.
1	May 15, 1998	wek	Updated lock out/tag out info
2	May 13, 1999	R. Hawthorne	Updated sec 2-3-3 DAT drive replacement with jumper and switch setting illustration.
3	May 21, 1999	S.M.Atladottir	Updated Procedure References for New GUI