

TABLE OF CONTENTS

TABLE OF CONTENTS	1
DESCRIPTION	2
1- ACCESSING THE INTERCOM BOARD	2
2- KEYBOARD SPEAKER VOLUME	3
2-1 Maximum Keyboard Speaker Volume R3 or R80 - (Patient's Voice).....	3
2-2 Minimum Keyboard Speaker Volume- R10 or R9 (Patient's Voice).....	5
3- PATIENT SPEAKER VOLUME (SPEAKER IN BORE)	5
3-1 Minimum Patient Speaker- Adjustment of R5 or R81 (Operator's voice).....	5
4- SYSTEM RESTORATION	6
A-APPENDIX.....	7
A-1 Intercom Board Version 2117167-4	7
REVISION HISTORY	8

DESCRIPTION

There are three volume adjustments for the Operator Workspace Keyboard Assembly: The intercom board provides additional volume adjustments:

- Keyboard speaker maximum volume
- Keyboard speaker minimum volume
- Patient speaker minimum volume

CAUTION

One of two models of the console intercom boards can be found in the system. 2117167-4 or 2204382. The model number can be found printed on the board itself. Insure you know what board you have before making adjustments.


CAUTION

One of two models of the remote intercom can be found in the system. 46-264758g1 or 2204384. The model number can be found printed on the board itself. Insure you know what board you have before making adjustments.

(These procedures require two people to perform them.)

1- ACCESSING THE INTERCOM BOARD

Access the Intercom Board per Illustration 1-1.

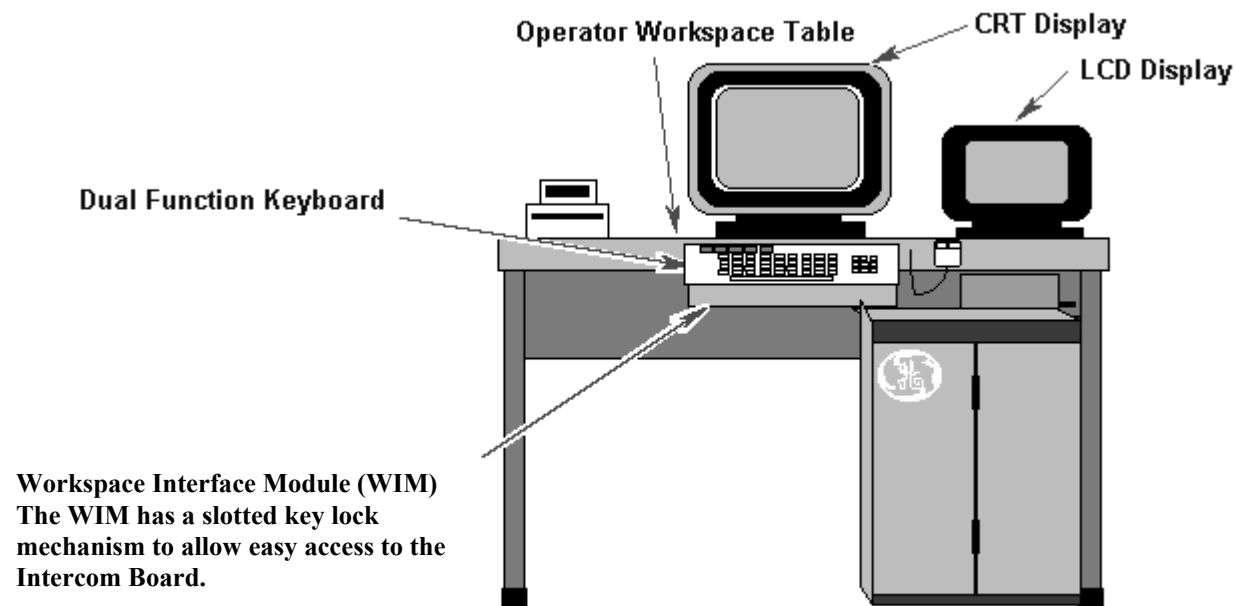

ACCESS THE INTERCOM BOARD
ILLUSTRATION 1-1

For instructions on removing the WIM from the workspace table and replacing it, see the document Workspace Table Replacements Octane, *Removal Procedures and Installation Procedures*.

2- KEYBOARD SPEAKER VOLUME


Note

There are jumpers on this board that can be moved to change the “default” amplifier output for both transmit (Microphone sensitivity) and Receive (Speaker volume) settings. It is recommended that the jumpers remain in the default position. However, if the customer wants more or less volume than the current volume levels provide, refer to Appendix A to adjust the Base Amplification of the board.


2-1 Maximum Keyboard Speaker Volume R3 or R80 - (Patient's Voice)

This section adjusts the maximum volume of the Patient's voice, in the magnet bore. The Operator benefits from this adjustment.


- Operator Volume- controls the speaker volume in the Bore, (Operator's voice).
- Patient Volume- controls the volume of the speaker in the keyboard, (Patient's voice). Refer to illustration 2-1.

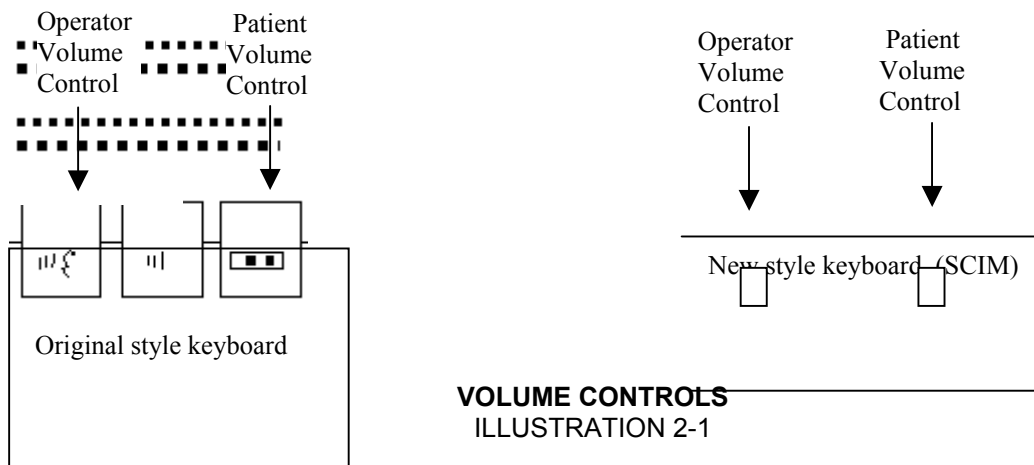

1. At the Operator Workspace keyboard, set the Patient Volume Control, (Middle control) at maximum, by rotating the volume wheel UP with your finger, (Clockwise).
2. One-person talks into the patient microphone in the magnet bore. The second person adjusts R3 or R80 on the Intercom Board for the desired maximum volume at the keyboard speaker. This is an iterative process. Repeat the voice test until the maximum desired level is achieved.
3. Depending on the version of Intercom board, adjust R3 or R80 clockwise to increase maximum volume at the patient speaker. Refer to illustrations 2-2 and 2-3 to determine which version of Intercom board you have.

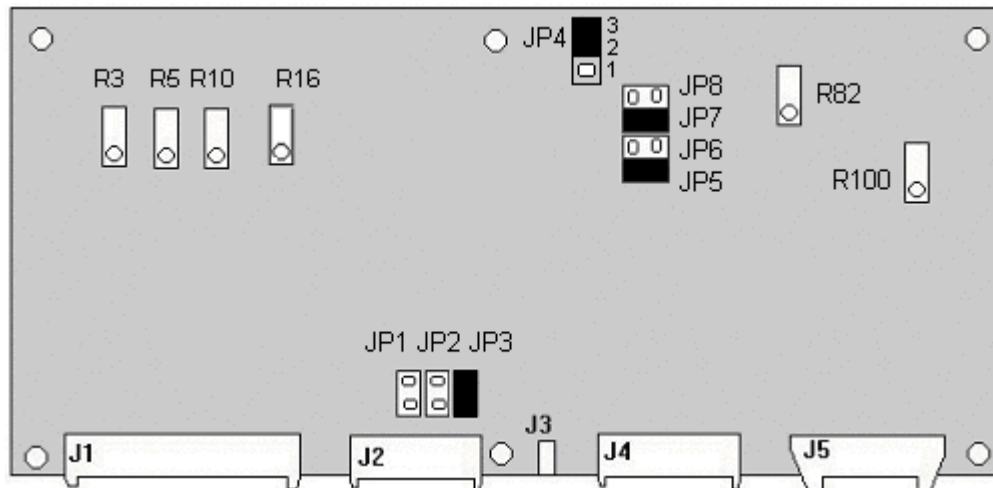

2117167-4 CONSOLE INTERCOM BOARD
ILLUSTRATION 2-2

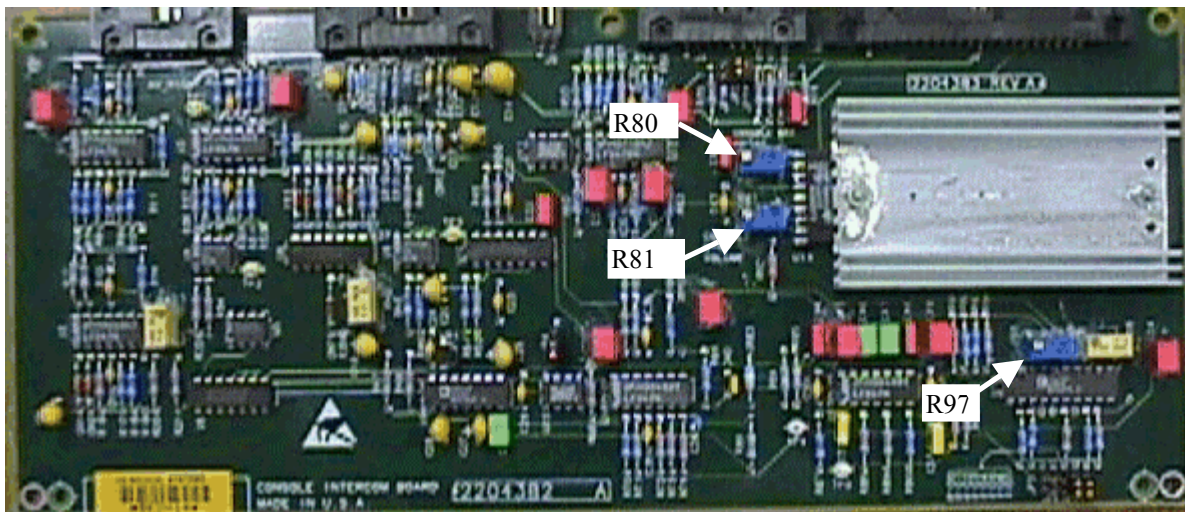

2204382 CONSOLE INTERCOM BOARD
ILLUSTRATION 2-3

Note

On the 2117167-4 console intercom board, R3 (keyboard speaker maximum volume) is factory set at approximately 11.5 turns clockwise. Turning R3 clockwise increases the maximum volume. (While turning the potentiometer, you will hear a faint clicking sound when it reaches the maximum or minimum limit.)

Note

On the 2204382 console intercom board, R80 (keyboard speaker maximum volume) is factory set at approximately 7.5 turns clockwise. Turning R80 clockwise increases the maximum volume. (While turning the potentiometer, you will hear a faint clicking sound when it reaches the maximum or minimum limit.)

2-2 Minimum Keyboard Speaker Volume- R10 or R9 (Patient's Voice)

This section adjusts the minimum volume of the Patient's voice, in the magnet bore. The Operator benefits from this adjustment.

Note

On the 2117167-4 Console Intercom Board, R10 (keyboard speaker minimum volume) is factory set at approximately ten turns clockwise.

- Turning R10 counterclockwise increases the minimum volume.

On the 2204382 Console Intercom Board, R97 (keyboard speaker minimum volume) is factory set at approximately 5.5 turns counter clockwise.

- Turning R97 counter clockwise increases the minimum volume.

While turning the potentiometer, you can hear a faint clicking sound when it reaches either the maximum or minimum limit.

1. At the Operator Workspace keyboard, set the Patient Volume Control, (Middle control) at minimum, by rotating the volume wheel down with your finger. (counterclockwise).
2. One person should listen at the keyboard speaker (Patient's Voice). The other person talks into the patient microphone in the magnet bore. The person at the keyboard should make minor adjustments to R10 or R97 on the Intercom Board for the desired minimum volume at the keyboard speaker. This is an iterative process. Make a slight adjustment to the potentiometer, test the volume and repeat the voice test until the desired Minimum Volume level has been achieved. Refer to Illustration 2-2 for location of R10 on board version 2117167-4, or Illustration 2-3 for location of R97 on board version 2204382.
3. Adjust R10 or R97 counterclockwise to increase minimum volume at the keyboard speaker.

3- PATIENT SPEAKER VOLUME (SPEAKER IN BORE)

This section adjusts the minimum volume of the Operator's voice, for the patient. The Patient lying in the magnet bore benefits from this adjustment.

3-1 Minimum Patient Speaker- Adjustment of R5 or R81 (Operator's voice)

1. At the Operator Workspace Keyboard, set the patient volume control at minimum by rotating the Operators volume wheel Down with your finger, (Counterclockwise). Refer to illustration 2-1. This control is to the far left.
2. One person listens at the patient speaker. The other person talks into the keyboard microphone and adjusts R5 or R81 on the Intercom Board for the desired minimum volume at the patient speaker. This is an iterative process and requiring repeated voice testing and potentiometer adjustment. Repeat the process until the desired minimum level is achieved.

Note

On the 46-264758G1 Remote Intercom Board, R5 (patient speaker minimum volume) is factory set at approximately four turns clockwise.

- Turning R5 clockwise increases the minimum volume.

On the 2204382 Console Intercom Board, R81 (patient speaker minimum volume) is factory set at approximately 7.5 turns clockwise. Turning R81 clockwise increase the minimum volume.

While turning the potentiometer, you can hear a faint clicking sound when it reaches the maximum or minimum limit.

4- SYSTEM RESTORATION

1. Close the Workspace Interface Module. Be sure that the module is completely sealed.
2. Perform a head and body scan to ensure system functionality.

A-APPENDIX

A-1 Intercom Board Version 2117167-4

Some customers may want more or less volume from the Intercom system. This appendix explains how to configure the jumpers on Intercom board version 2117167-4, to increase or decrease Base Amplification of the board.

Note

There are jumpers on this board that can be moved to change the “default” amplifier output for both transmit (Microphone sensitivity) and Receive (Speaker volume) settings. It is recommended that the jumpers remain in the default position. However, if the customer wants more or less volume than the current volume levels provide, reset the jumpers accordingly, and then proceed to fine tuning using the potentiometers as described in the next few paragraphs. Refer to Table A-1.

TABLE A-1
DEFAULT JUMPER SETTINGS
2117167-4

Jumper	JP1	JP2	JP3	JP4	JP5	JP6	JP7	JP8
position	Off	Off	On	On	On	Off	On	Off
Pin #				2-3				

Base Amplification Choices:

- Jumpers JP5 and JP6 adjust the Base Amplification for the Patient Speaker located in the KEYBOARD. If the Operator wishes to hear the patient better, move the Jumper on JP5 to JP6, and readjust the potentiometers for Minimum and Maximum volume levels.
- Jumpers JP7 and JP8 adjust the Base Amplification for the Operators Speaker, located within the MAGNET BORE. If the Patients are complaining that they cannot hear the Operator, move the Jumper on JP7 to JP8, and readjust the potentiometers for Minimum and Maximum volume levels.

Proceed with adjustments as described in section 2-1 through 3-1.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
A	Apr 3, 1998	R. Hawthorne	Initial conversion from toolbook. Added New adjustments for 2204382 console interface board.
0	Aug 7, 1998	M. Keber	Corrections to locations of speaker and listener for Keyboard speaker volume adjustments. Miscellaneous style guide cleanup for procedure.
1	Dec 12, 2002	D. Hofstetter	Cleared up adjustment process and added graphics and appendix.