

## TABLE OF CONTENTS

|                                                            |   |
|------------------------------------------------------------|---|
| TABLE OF CONTENTS .....                                    | 1 |
| DESCRIPTION .....                                          | 2 |
| 1- INTRODUCTION .....                                      | 2 |
| 2- TEST SET-UP .....                                       | 2 |
| 2-1 Remove Power From The Operator Workspace Cabinet ..... | 2 |
| 2-2 Workspace Part Removal.....                            | 2 |
| 3- VOLTAGE ADJUSTMENTS.....                                | 4 |
| 4- PUTTING SYSTEM BACK TOGETHER .....                      | 5 |
| REVISION HISTORY.....                                      | 6 |

## DESCRIPTION

This procedure describes how to adjust power supply voltages for the Workspace Interface Module (WIM) and the SCSI Expansion Box (PED) located in the AC Power Distribution Box in the computer cabinet.

## 1- INTRODUCTION

The voltages are pre-adjusted when the workspace cabinet is assembled. All replacement parts will have the voltage adjusted to the proper level. Therefore, this procedure will not be used often.



**The voltage adjustments must be made with power on and the WIM disassembled.**

## 2- TEST SET-UP

### 2-1 Remove Power From The Operator Workspace Cabinet




**FATAL ELECTRIC SHOCK HAZARD!! ENSURE THAT POWER IS OFF TO THE OPERATOR WORKSPACE BEFORE STARTING ANY PROCEDURE IN THIS DOCUMENT.**

1. Log out from Host Computer.
2. Shut down the Host Computer.
3. Select **[Start]** from the lower left corner and select "Shutdown" from the pull-down menu. Then, select "Shutdown the Computer" and press **[Yes]**.
4. As applicable, open the circuit breaker(s) that power the Operator Workspace to turn power off.
5. Lock and tag the circuit breaker(s).

### 2-2 Workspace Part Removal

1. Remove the SCSI Expansion Box. Refer to the procedure Workspace Cabinet Replacements Octane. The adjustment points will be exposed on the AC Distribution Box.
2. Carefully re-apply power to the Workspace Cabinet.
3. Turn off power to the Workspace Interface Module and SCSI Expansion Box by pressing the power buttons on the front of the Workspace cabinet. Refer to Illustration 2-1.

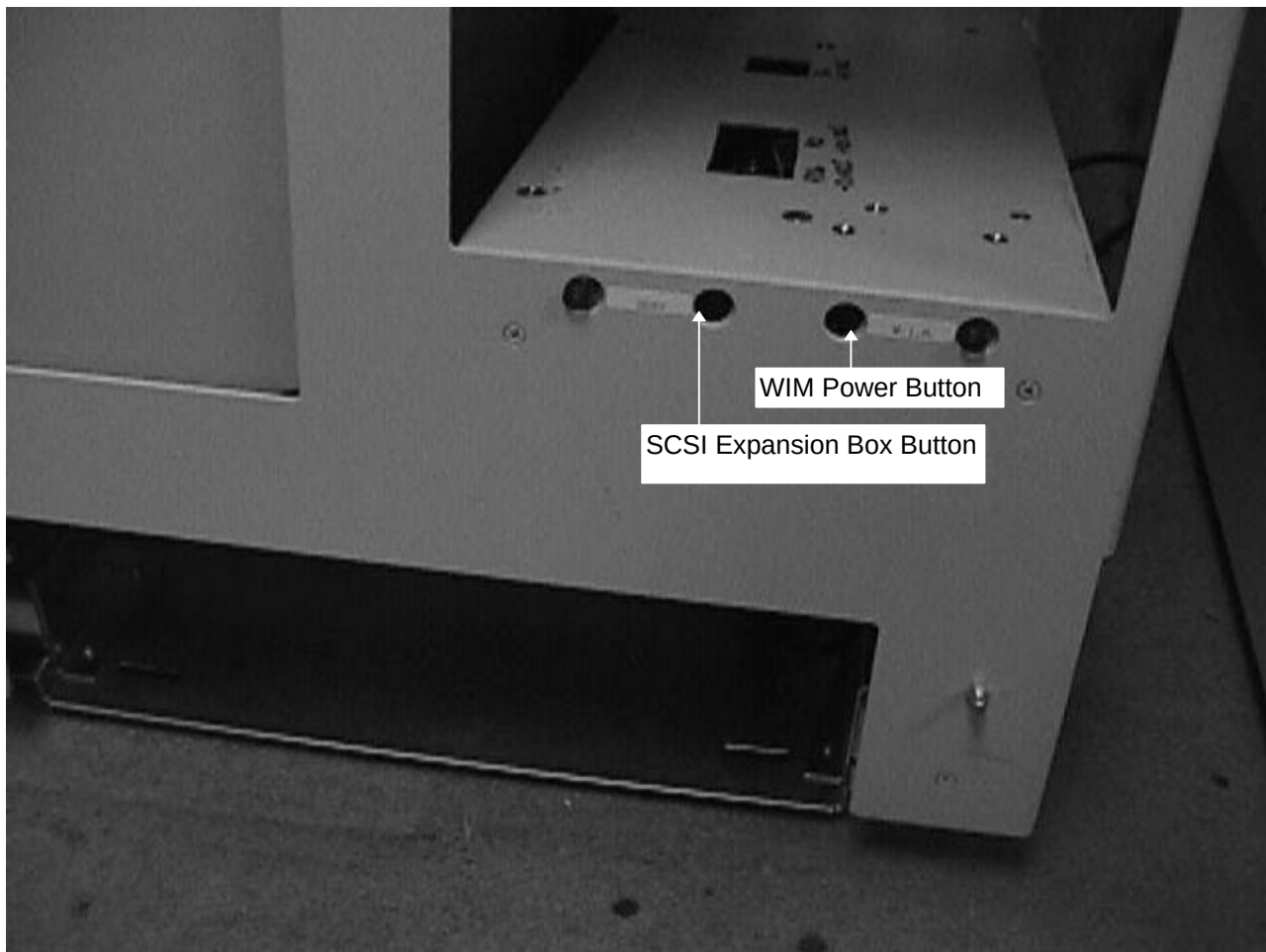

**WIM POWER BUTTON**  
ILLUSTRATION 2-1

4. Remove the keyboard cable from the Workspace Interface Module (WIM) at the back of the table.
5. Slide out the WIM board using the handle releases and place on convenient location to connect to test points.



**Circuit boards are static sensitive. The static sensitive components on these circuit boards may be damaged if not handled in a static free environment. Take appropriate precautions (e.g. wear properly grounded wrist strap) when handling circuit boards.**


6. Remove the Intercom board (top board) on the WIM board.
7. Carefully, re-apply power to the WIM by pressing the power switch on the front of the Workspace Cabinet.

### 3- VOLTAGE ADJUSTMENTS

1. The voltage adjustments are done via the potentiometers in the AC Distribution Box. Refer to Illustration 3-1.

**Note**

The potentiometer labeling may appear differently on your box.

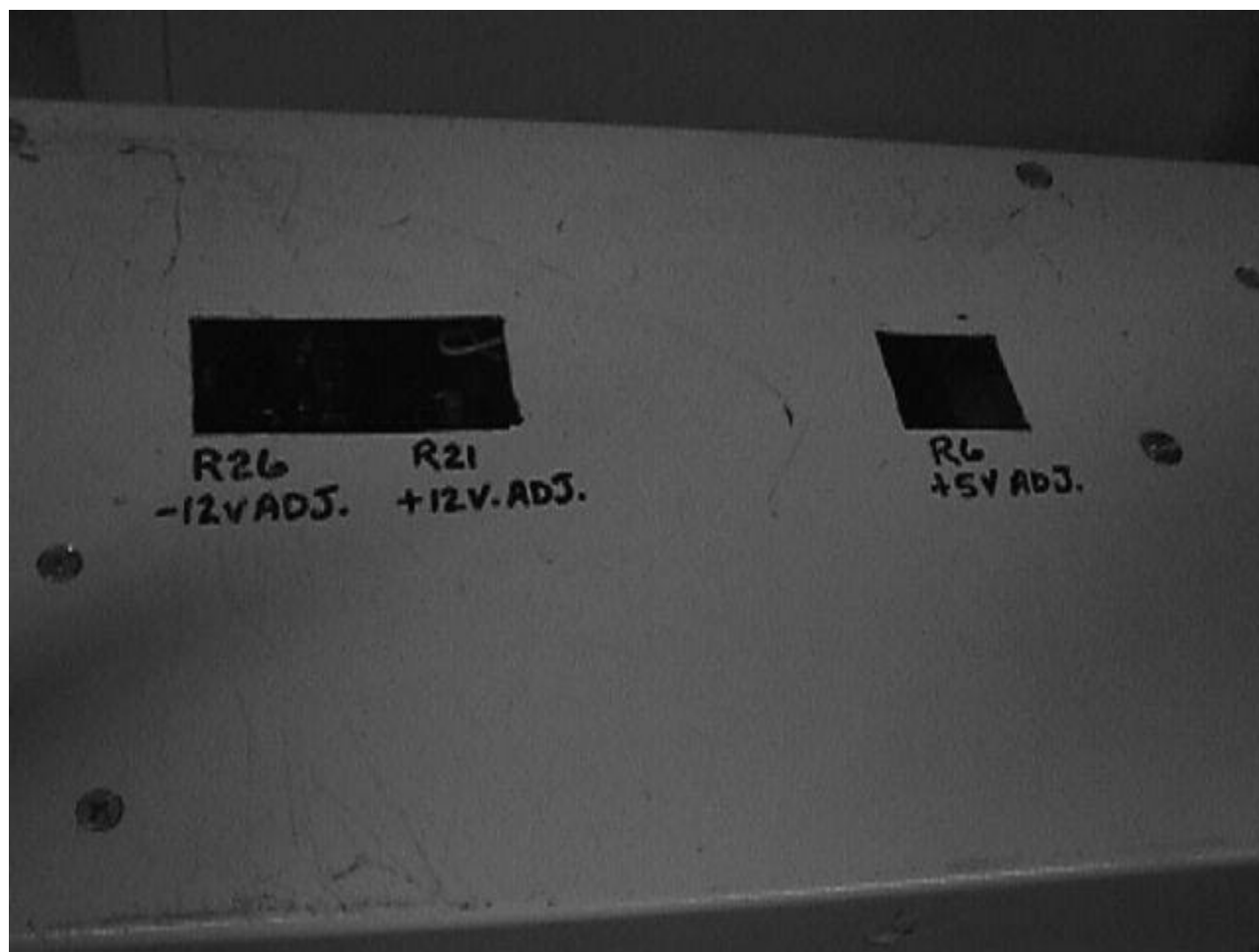

**VOLTAGE ADJUSTMENT POTS**  
ILLUSTRATION 3-1

2. The test points are located on the WIM. Use Table 3-1 for the test points and voltage adjustment levels.

**TABLE 3-1**  
**VOLTAGE ADJUSTMENTS**

| Voltage | Potentiometer | Test Point (+) | Test Point (Common) | Adjustment Level  |
|---------|---------------|----------------|---------------------|-------------------|
| +12 VDC | R26           | TP2            | TP1 or TP37         | 12.0 +/- .02 VDC  |
| -12 VDC | R21           | TP3            | TP1 or TP37         | -12.0 +/- .02 VDC |
| +5 VDC  | R6            | TP29           | TP28                | 5.10 to 5.12 VDC  |

**4- PUTTING SYSTEM BACK TOGETHER**

1. Remove power from the WIM and SCSI Expansion Box by pressing the power switch on the front of the Workspace Cabinet.
2. Remove power to the Workspace Cabinet.
3. Re-connect the top board on the WIM and install the WIM back inside the module on the Workspace Table.
4. Re-connect the keyboard cable to the WIM.
5. Return the SCSI Expansion Box to its proper location and re-connect all cables.
6. Return power back to the system.
7. Perform at least one test scan to verify the system has been re-connected properly

## REVISION HISTORY

| REV | DATE        | AUTHOR | PRIMARY REASONS FOR CHANGE |
|-----|-------------|--------|----------------------------|
| 0   | May 4, 1998 | K. L-P | Initial revision           |
|     |             |        |                            |
|     |             |        |                            |
|     |             |        |                            |
|     |             |        |                            |