

TABLE OF CONTENTS

table of contents.....	1
1- Tools Required	2
2- Removal of Failed Coupler.....	2
3- Installation of New Coupler.....	3
4- Calibration Requirements	3
5- Functional Check Requirements	3
revision history.....	4

Description - This section applies to Signa systems installed with the 1.0T and 1.5T RF/Pen II cabinet.

1- TOOLS REQUIRED


#2 Phillips screwdriver

WARNING!

POSSIBLE PERSONAL INJURY! REMOVE POWER FROM THE CABINET BEFORE REMOVING ANY COMPONENTS FROM THE CABINET. VERIFY THAT ALL BREAKERS IN THE REAR OF CABINET ARE OFF, AND DISCONNECT MAIN POWER FEED CORD BEFORE OPENING CHASSIS FOR SERVICE. SERIOUS INJURY OR DEATH BY ELECTROCUTION MAY OTHERWISE OCCUR. LOCK OUT AND TAG OUT THE CABINET.

2- REMOVAL OF FAILED COUPLER

1. Turn off the cabinet and RF amplifier circuit breakers on the rear of the cabinet. Lock out and tag out the cabinet (see Illustration 2-1). (Refer to *Procedure For Safety: Section 6.*)

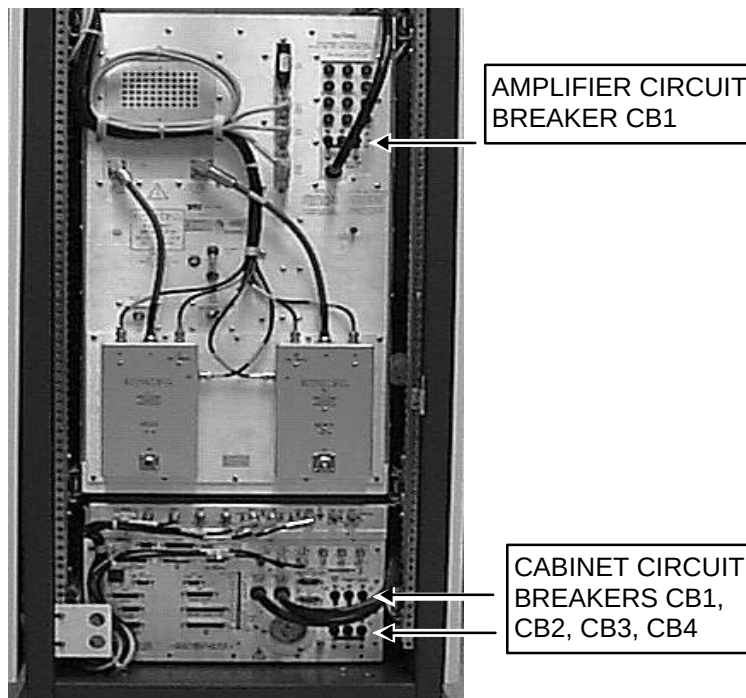


LOCK OUT AND TAG OUT
ILLUSTRATION 2-1

2. Verify that power to the RF/Pen II Cabinet is off by checking the fans and LEDs on both the RF amplifier and the System Support Module (SSM).
3. Disconnect the eight cables attached to the quad directional coupler.

2- REMOVAL OF FAILED COUPLER (continued)

4. Disconnect the two cables that attach the quad directional coupler to the RF amplifier.

5. Remove the ten 8-32x7/16 screws that secure the quad directional coupler to the RF amplifier (see Illustration 2-2).

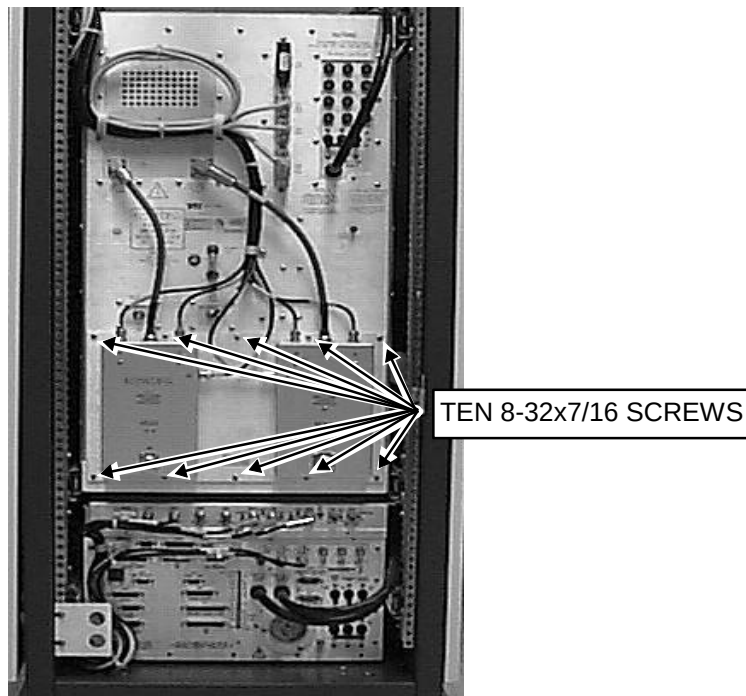

REMOVAL OF QUAD DIRECTIONAL COUPLER
ILLUSTRATION 2-2

6. Remove the quad directional coupler from the rear of the cabinet.

3- INSTALLATION OF NEW COUPLER

1. Position the new quad directional coupler on the rear of the cabinet.
2. Secure, using ten 8-32x7/16 screws.
3. Connect the eight cables to the quad directional coupler, matching the labels on the cables to the connections.
4. Connect the two cables of the quad directional coupler to the RF amplifier.
5. Remove lock out and tag out, and restore power to the cabinet and RF amplifier.

4- CALIBRATION REQUIREMENTS

Perform APM calibration (refer to the procedure for *APM Calibration*).

5- FUNCTIONAL CHECK REQUIREMENTS

Perform one head and body scan.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	June 15, 1998	Erbtec Engineering	Initial Toolbook conversion to Word. (KK)
1	MAY 21, 1999	SMA	Updated Procedure References for New GUI