

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1 - OVERVIEW AND FLOWCHART	2
2 - BODY AND HEAD MAXIMUM POWER RF OUTPUT CHECK	5
2-1 Tools and Instruments Required	6
2-2 System Preparation	7
2-3 Body RF Output Power Check	8
2-4 Head RF Output Power Check	12
3 - BODY AND HEAD RF OUTPUT POWER ADJUSTMENTS	15
3-1 Body Gain Pot Adjustment	15
3-2 Head Gain Pot Adjustment	20
4 - SYSTEM RESTORATION	24
APPENDIX A — ALTERNATE EQUIPMENT SETUPS & NON-SERVICE PROTOCOLS	25
APPENDIX B — CALIBRATION REQUIREMENTS AFTER HARDWARE REPLACEMENT	28
APPENDIX C — SOFTWARE AND FORMULAE FOR MISCELLANIOUS RF CONVERSIONS	29
C-1 RF Conversion Software	29
C-1-1 RF Power Calculator	29
C-1-2 Power Calculator Tool	29
C-2 RF Conversion Formulae	29
C-3 1.5T Division Conversions	30
APPENDIX D — CHARACTERIZATION FOR SCOPES WITH < 300 MHZ BANDWIDTH	31
D-1 Scope Channel Characterization Done During Calibration	31
D-2 Scope Channel Characterization Done By Comparison	31
D-3 Scope Characterization Using the RF Power Measurement Kit Signal Reference	31
APPENDIX E — (U)CERD EXCITER OR RRF MASTER EXCITER RF SIGNAL OUTPUT	33
E-1 Exciter or RRF Master Exciter RF Output Measurement Using an Oscilloscope	34
E-2 Exciter or RRF Master Exciter RF Output Check Using The RF Power Measurement Kit	37
APPENDIX F — DUMMY LOAD AND CABLES CALIBRATION	40
F-1 Overview	40
F-2 Tools and Instruments Required	41
F-3 Initial Setup	41
F-4 Data Collection	44
F-5 Analysis	47
F-6 Calculation of RF power	48
F-7 System Restoration	49
F-8 Wattmeter Trivia	50
F-9 Scope Trivia	51
APPENDIX G — GENERAL TROUBLESHOOTING	53
REVISION HISTORY	56

1 - OVERVIEW AND FLOWCHART


This procedure describes how to calibrate the SRFD2 in a 1.5T system. Please see illustration 1-1 on the next page to confirm that the 1.5T SRFD2 cabinet pictured matches the cabinet to be checked and/or adjusted. The SRFD2 Cabinet contains the SRFD2 module and SSM but does not contain gradient hardware or a PDU. The SRFD2 should come from the factory pre-calibrated and, ideally, no initial calibration of the module in the field should be needed, however, this procedure must be performed to confirm that it is calibrated.

Calibration is done to prevent faults and to make sure that the SRFD2 is outputting the specified 16kW of RF body power and 2kW of RF head power. Adjustments may need to be made if any of the hardware in the RF chain has been replaced (See Appendix B), output levels have decreased with age, or if any of the measurements taken during the check are found not to meet specifications. The head and body RF power levels will be checked in section 2 and, if in specification, no adjustments will be made. The head and body power levels and the gain adjustments are made in section 3. Please see the flowchart in Illustration 1-2. Appendices A - G are referred to in the procedure and provide extra information and directions to assist the user with accomplishing the RF power measurement task. Appendix G provides helpful troubleshooting information in the event that a problem is encountered during the RF power check or calibration.

Note


The Body and Head Gain pots are in series with each other. The Body Gain pot is first in the series. The Head Gain pot is factory set so that the head gain is 9 dB less (63 dBm) than the body gain (72 dBm).

1 - OVERVIEW AND FLOWCHART (CONTINUED)

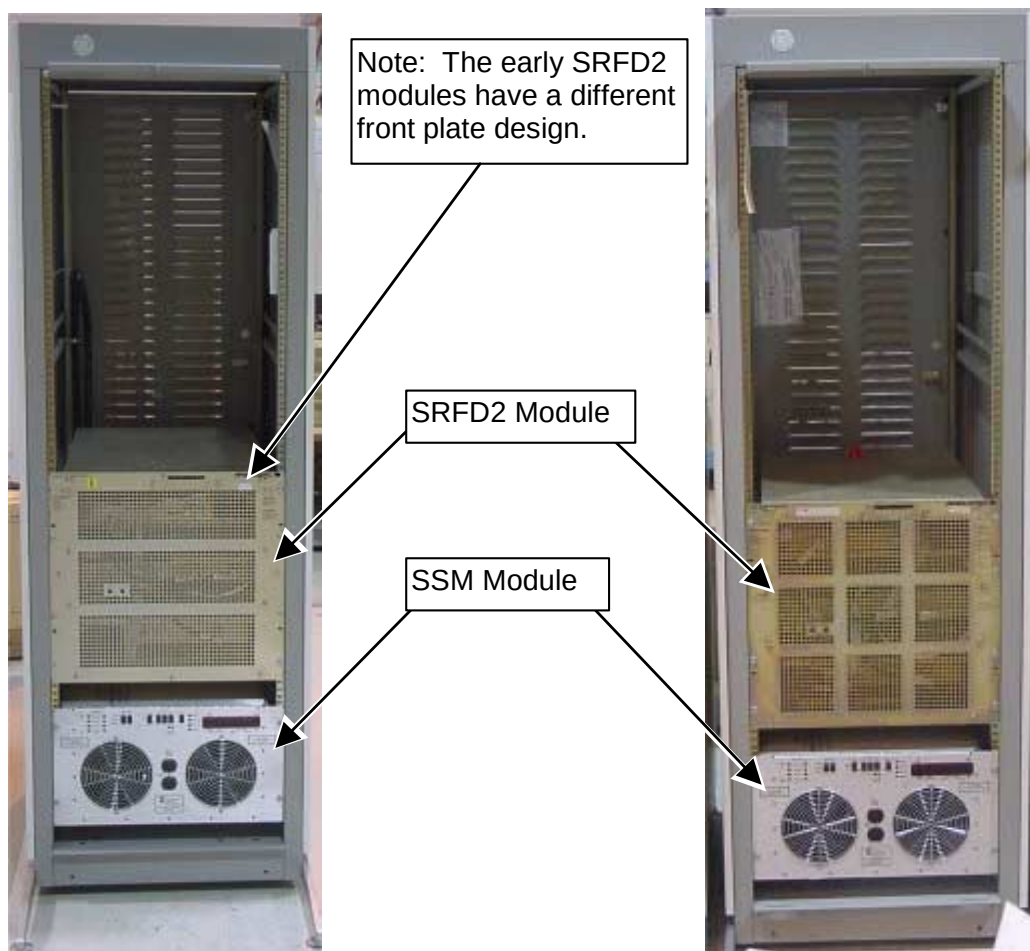

1.5T SRFD2 CABINET
ILLUSTRATION 1-1

1 - OVERVIEW AND FLOWCHART (CONTINUED)

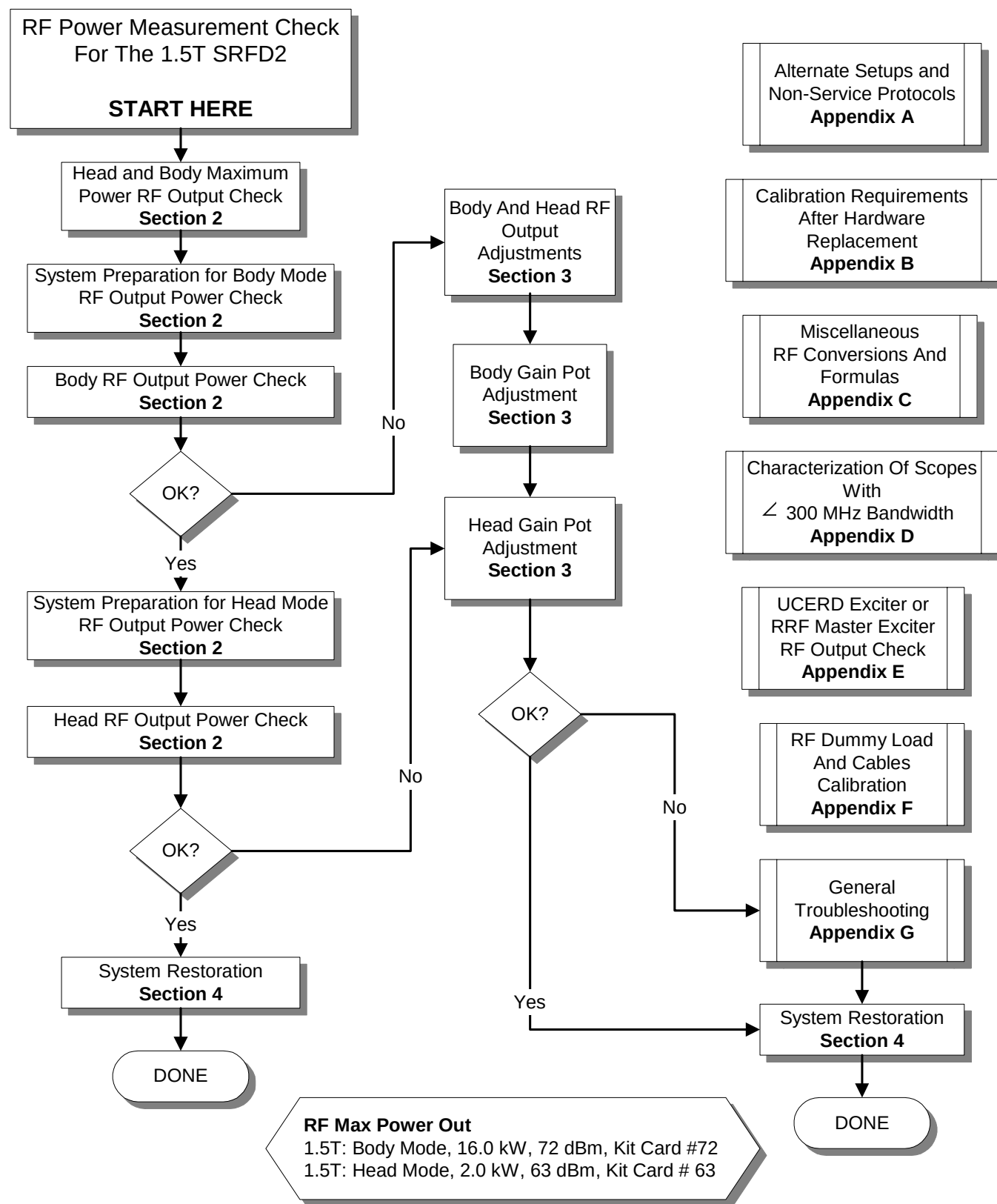

CALIBRATION FLOWCHART
ILLUSTRATION 1-2

2 - BODY AND HEAD MAXIMUM POWER RF OUTPUT CHECK

This process will test the assumptions that:

- The (U)CERD Exciter or RRF (Remote RF) Master Exciter is outputting enough power for the system to meet the rated head and body RF output power levels (Appendix E describes how to check the (U)CERD Exciter or RRF Master Exciter output.).
- The SRFD2 has been pre-adjusted at the factory and that, ideally, no further adjustments by the FE are necessary.
- The RF amplifiers within the SRFD2 are capable together of outputting 16kW.

The (U)CERD Exciter or RRF Master Exciter cannot be adjusted. Appendix E describes how to measure and verify that the Exciter or RRF Master Exciter output is correct. The head and body RF output power levels must be checked and, if not correct, adjusted into specification. Three methods are available and listed in the order of preference for measuring head and body RF output power:

- **RF Power Measurement Kit** - Easy to use, preferred method of power measurement.
- **Bird wattmeter** - Must account for cable and load loss to determine actual power.
- **Oscilloscope** - Must know the scope channel correction factor if the scope bandwidth is < 300 MHz. Must account for cable and dummy load loss and the scope bandwidth limitation (channel correction factor) to accurately calculate the power from the observed peak voltage. This method works but the process used to derive the power is somewhat complex and lacks the accuracy of the previous two methods. It is, for this reason, not recommended.

It is *strongly recommended* that the RF Power Measurement Kit be used to obtain an accurate measurement of the RF output power. If it is not possible to obtain an RF Power Measurement Kit then one of the two other methods can be used provided that:

- The Dummy Load and Cables Calibration procedure in Appendix F has been performed and the individual loss values are known.
- The oscilloscope input channel correction factor is known when the scope bandwidth is < 300 MHz *and* the wattmeter is not being used. See Appendix D.

WARNING!

FAILURE TO ACCOUNT FOR THE ACTUAL LOSSES IN THE DUMMY LOAD, CABLES, AND OSCILLOSCOPE WHEN NOT USING THE RF POWER MEASUREMENT KIT TO MEASURE RF POWER CAN RESULT IN SIGNIFICANT MEASUREMENT ERRORS.

A Calibration Flowchart, Illustration 1-2, has been provided in the Overview, Section 1.

2-1 Tools and Instruments Required

Refer to Table 2-1 for equipment required to measure power with the RF Power Measurement Kit or refer to Appendix A (Alternate Equipment Setup). The RF Power Measurement Kit is the preferred method for RF power measurement.

TABLE 2-1
EQUIPMENT REQUIRED IN ADDITION TO THE RF POWER MEASUREMENT KIT


Item	Description	Part Number
1.	100 MHz Scope (equivalent or greater)	46-183029P61
2.	RF Power Measurement Kit	46-317724G1 or G2
3.	50 ohm, 200 Watt, 30dB Attenuator (dummy load)	46-317724P14
4.	Screwdriver, slotted (long & narrow)	—
5.	RF Test Cables Kit	46-255816G1
6.	Digital Multimeter (DMM)	46-194427P49
7.	DMM Test Cable – BNC to Banana (optional)	2239139

2-2 System Preparation



PROPERTY DAMAGE! PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, BY REMOVING ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

1. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the **DANGER** message on this page.
2. Remove the front cabinet cover and open the rear cabinet cover.
3. See Illustration 2-1. On the front of the SSM place the:
 - 2 (two) power MONITOR switches (A and B) to the middle BYPASS position.
 - 3 (three) DRV FLT switches to the bottom DIS (disable faults) position.

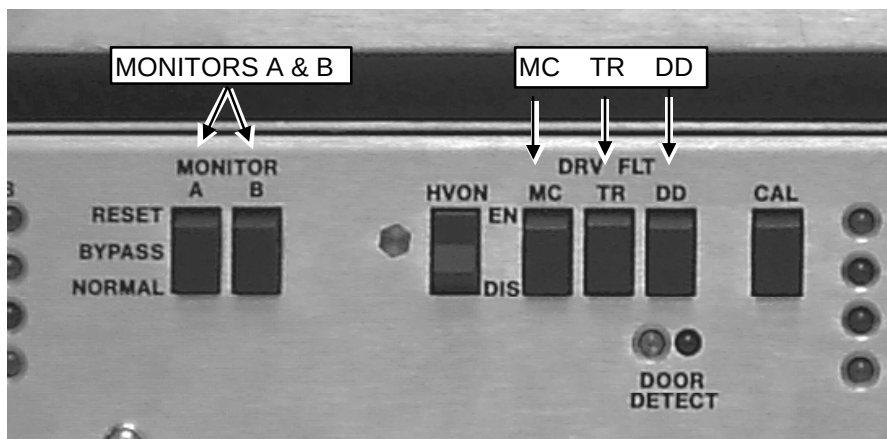

SSM FRONT PANEL SWITCHES DISABLED
ILLUSTRATION 2-1

4. **If you are using the RF Power Measurement Kit** then calibrate the scope by referring to the RF Power Measurement Kit laminated card set.
 - a. Look in the upper right corner of each card and find the card labeled **CAL**.
 - b. Configure the scope as in the illustration on the card.
 - c. Follow the directions on the card to calibrate the scope using the 4 dBm calibrator.
5. **If you are not using the RF Power Measurement Kit** then complete the Dummy Load and Cables Calibration procedure in Appendix F. If using the scope to measure the RF output power then also make sure that the scope correction factor has been determined for the input channel to be used as described in Appendix D

2-3 Body RF Output Power Check



PROPERTY DAMAGE! PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, BY REMOVING ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

1. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the two **DANGER** messages on this page.



PERSONAL INJURY! PREVENT POSSIBLE RF BURNS WHEN DISCONNECTING HELIAX CABLES FROM J3 OR J4 ON THE REAR OF THE SRFD2 BY VERIFYING THAT THE SYSTEM IS NOT MANUALLY PRESCANNING OR SCANNING. VERIFY THAT THE SCAN DESKTOP ICON DISPLAYS THE "IDLE" MESSAGE.



2-3 Body RF Output Power Check (continued)

2. **If using the RF Power Measurement Kit** then refer to the RF Power Measurement Kit laminated card set.
 - a. Look in the upper, right corner of each card and find the card labeled **72** (1.5T Body Output).



The body RF output is no longer connected to the non-existent EFB unit, as the older RF Power Measurement Kit cards show, but instead to the J4 output on the rear of the SRFD2 module. An RF adapter is provided in the RF Power Measurement Kit to connect between the HN J4 body RF output and the RF Power Measurement Kit 40 dB N-connector coupler.

- b. Configure the system as shown in the illustration on the card.
 - c. Confirm that the rotary attenuator is set to the correct position indicated on the card.
3. **If using the wattmeter or scope** (NOT the RF Power Measurement Kit) to measure power then refer to Appendix A (Alternate Equipment Setup) for the proper system body configuration.

2-3 Body RF Output Power Check (continued)

4. Prepare the system to scan in Body mode per Table 2-3 or refer to Appendix A (Non-Proprietary protocol).

Note

If the system is operating with Release 9.X, CNV4 software, or its equivalent then it is no longer necessary to modify the **ia_rf1** and **ia_rf2** CV values. These will be set automatically when the **calmode** CV value is set equal to 5.

TABLE 2-3
SCAN PROTOCOL: BODY MODE

Note: This is the alternate proprietary procedure available for GE use and for sites with a valid Advanced Service Package Limited License.
Refer to Appendix A for the non-proprietary protocol.

A. **[New Pt]**

Id: **geservice<ENTER>**

Name: **rfi cals**

Weight (Lb.): **300<ENTER>**

Set Patient Protocols to **Service** from the pull-down menu.

B. At front enclosure:

Landmark in the Head area—remove any coils.

press **LANDMARK**.

press **MOVE TO SCAN**.

C. In the Patient Position Protocol field:

type **o.41.1<ENTER>** (o=Other, 41.1 =series) to load the body protocol

OR select **other** and select protocol **41** and select series **1**.

D. **[Save Series]**.

E. **[Research Operations]**.

[Setup Params]. Set TG to **50 [Done]**.

F. **[Research Operations]**.

[Display CVs]. Highlight CV Name and enter the following:

CV Name: **calmode<ENTER>**, CV Value: **5<ENTER>** (Dual Logamp Waveform).

Note

Skip the next two steps if the system has Release 9.X, CNV4, or equivalent software. The software will set the **ia_rf1** and **ia_rf2** CV values automatically once the **calmode** CV Value is set equal to 5.

CV Name: **ia_rf1<ENTER>**, CV Value: **32766<ENTER>** (sets 90° pulse full scale).

CV Name: **ia_rf2<ENTER>**, CV Value: **0<ENTER>** (turns off 180° pulse).

[Accept].

G. **[Research Operations]**

H. **[Download]**.

5. Verify Body LED is illuminated on front of SRFD2 Module.
6. **[Manual Prescan] [Scan TR]**. View the Body RF output waveform on the oscilloscope. Slowly increase TG to 200.

2-3 Body RF Output Power Check (continued)



Check the system center frequency NOW, especially if this is a new installation, and confirm that it consists of 8 digits and that it is correct. The RF amplifiers will output a low and distorted RF waveform if center frequency is wrong due to a missing digit.

7. Verify the Unblank LED on the front of the SRFD2 Module pulses on/off indicating that the amplifier is outputting RF energy.
8. If the SRFD2 faults with error 2254212 (peak power fault) and the cabling to and from the SRFD2 is correct, then it will be necessary to calibrate the amplifier. Proceed to **Section 3 - Body And Head RF Output Power Adjustments**.
9. **If using the RF Power Measurement Kit** then refer to the specified number of divisions needed to achieve 16kW output from the RF Power Measurement Kit card **72** (1.5T Body RF Output) and verify that the waveform displayed on the scope meets, and does not exceed, this specification.
10. **If using the wattmeter procedure** then read the wattmeter display and use the formula below to calculate the body RF power and verify that it meets, and does not exceed, the 16kW (72dBm) specification. The dummy load and cable loss factor was determined from the procedure in Appendix F.

RF Power Measurement (in watts) Using Wattmeter And Formula:

Wattmeter reading (in watts) **X** dummy load and cable loss factor

11. **If using the oscilloscope procedure (NOT the RF Power Measurement Kit)** then read the peak voltage (V_{peak}) from the scope display and use the formula below or the [Power Calculator](E:\rf\power\pwrcalc.htm) Tool (located at E:\rf\power\pwrcalc.htm on the Service Methods CD-ROM) to calculate the RF output power and verify that it meets, and does not exceed, the 16kW (72dBm) specification. The dummy load and cable loss factors were determined from the procedure in Appendix F. The scope correction factor was determined in Appendix D.

RF Power Measurement (in watts) Using Oscilloscope And Formula:

$$\frac{\left(\frac{V_{\text{peak}}}{\text{scope correction factor}} \right)^2}{100} \times \text{dummy load and cable loss factor}$$

12. Decrease TG to 0 (zero). **[Done] [End Exam]**.

2-3 Body RF Output Power Check (continued)

13. Verify that the scan desktop icon is displaying the "Idle" message.
14. If the body RF power output did not meet specifications then proceed directly to **Section 3 - Body And Head RF Output Power Adjustments**.
15. If the body RF power output met specifications then proceed to the next section, **Section 2-4 - Head RF Output Power Check**.

2-4 Head RF Output Power Check



PROPERTY DAMAGE! PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, BY REMOVING ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

1. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the two **DANGER** messages on this page.



PERSONAL INJURY! PREVENT POSSIBLE RF BURNS WHEN DISCONNECTING HELIAX CABLES FROM J3 OR J4 ON THE REAR OF THE SRFD2 MODULE BY VERIFYING THAT THE SYSTEM IS NOT MANUALLY PRESCANNING OR SCANNING. VERIFY THAT THE SCAN DESKTOP ICON DISPLAYS THE "IDLE" MESSAGE.



2. **If using the RF Power Measurement Kit** then refer to the RF Power Measurement Kit laminated card set.
 - a. Look in the upper, right corner of each card and find the card labeled **63** (1.5T Head Output).



The head RF output connection is no longer to the non-existent EFB unit, as the reference cards in some of the older kits show, but instead to the J3 output on the rear of the SRFD2 Module.

- b. Configure the system as shown in the illustration on the card.
 - c. Confirm that the rotary attenuator is set to the correct position indicated on the card.

2-4 Head RF Output Power Check (Continued)

3. **If using the wattmeter or scope** (NOT the RF Power Measurement Kit) to measure power then refer to Appendix A (Alternate Equipment Setup) for the proper system head configuration.
4. Prepare the system to scan in Head mode per Table 2-4 or refer to Appendix A (Non-Proprietary protocol).

Note

If the system is operating with Release 9.X, CNV4 software, or its equivalent then it is no longer necessary to modify the **ia_rf1** and **ia_rf2** CV values. These will be set automatically when the **calmode** CV value is set equal to 5.

TABLE 2-4
SCAN PROTOCOL: HEAD MODE

Note: This is the alternate proprietary procedure available for GE use and for sites with a valid Advanced Service Package Limited License.
Refer to Appendix A for the non-proprietary protocol.

- A. **[New Series]**
At Patient Protocols – select **other**.
- B. In the protocol field, type **o.41.2<ENTER>** (o=Other, 41.2 =series) to load the head protocol
OR select **[o.41] [Series 2] [Accept]**.
- C. **[OK]** (if required).
- D. **[Save Series]**
- E. **[Prepare to Scan]**
- F. **[Research Operations]**.
[Setup Params]. Set TG to **50 [Done]**.
- G. **[Research Operations]**.
[Display CVs]. Highlight CV Name and enter the following:
CV name: **calmode<ENTER>, 5<ENTER>** (Dual Logamp Waveform).

Note

Skip the next two steps if the system has Release 9.X, CNV4, or equivalent software. The software will set the **ia_rf1** and **ia_rf2** CV values automatically once the **calmode** CV Value is set equal to 5.

CV name: **ia_rf1<ENTER>, 32766<ENTER>** (sets 90° pulse full scale).

CV name: **ia_rf2<ENTER>, 0<ENTER>** (turns off 180° pulse).

[Accept].

- H. **[Research Operations]**
- I. **[Download]**.

5. Verify Head LED is illuminated on front of SRFD2 Module.
6. **[Manual Prescan] [Scan TR]**. View the Head RF output waveform on the oscilloscope. Slowly increase TG to 200.

2-4 Head RF Output Power Check (Continued)



Check the system center frequency **NOW**, especially if this is a new installation, and confirm that it consists of 8 digits and that it is correct. The RF amplifier will output a low and distorted RF waveform if center frequency is wrong due to a missing digit.

7. Verify the Unblank LED on the front of the SRFD2 Module pulses on/off indicating that the amplifier is outputting RF energy.
8. **If using the RF Power Measurement Kit** then refer to the specified number of divisions needed to achieve 2kW (63dBm) output from the RF Power Measurement Kit card **63** (1.5T Head RF Output) and verify that the waveform displayed on the scope meets, and does not exceed, this specification.
9. **If using the wattmeter procedure** then read the wattmeter display and use the formula below to calculate the head RF power and verify that it meets, and does not exceed, the 2kW (63dBm) specification. The cable loss factor was determined from the procedure in Appendix F.

RF Power Measurement (in watts) Using Wattmeter And Formula:

Wattmeter reading (in watts) **X** cable loss factor

10. **If using the oscilloscope procedure (NOT the RF Power Measurement Kit)** then read the peak voltage (V_{peak}) from the scope display and use the formula below or the [Power Calculator](E:\rf\power\pwrcalc.htm) Tool (located at E:\rf\power\pwrcalc.htm on the Service Methods CD-ROM) to calculate the RF output power and verify that it meets, and does not exceed, the 2kW (63dBm) specification. The dummy load and cable loss factor was determined from the procedure in Appendix F. The scope correction factor was determined in Appendix D.

RF Power Measurement (in watts) Using Oscilloscope And Formula:

$$\frac{\left(\frac{V_{\text{peak}}}{\text{scope correction factor}} \right)^2}{100} \text{ X dummy load and cable loss factor}$$

11. Decrease TG to 0 (zero). **[Done] [End Exam]**.
12. Verify that the scan desktop icon is displaying the "Idle" message.
13. If the amplifier is not outputting 2kW of RF power but the body power is within specifications then proceed directly to **Section 3-2 - Head Gain Pot Adjustment**.
14. If the RF head output power is 2kW and the body output power is 16kW then no further adjustments are necessary. Proceed directly to **Section 4 – System Restoration**.

3 - BODY AND HEAD RF OUTPUT POWER ADJUSTMENTS

3-1 Body Gain Pot Adjustment



PROPERTY DAMAGE! PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, BY REMOVING ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

1. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the two **DANGER** messages on this page.



PERSONAL INJURY! PREVENT POSSIBLE RF BURNS WHEN DISCONNECTING HELIAX CABLES FROM J3 OR J4 ON THE REAR OF THE SRFD2 MODULE BY VERIFYING THAT THE SYSTEM IS NOT MANUALLY PRESCANNING OR SCANNING. VERIFY THAT THE SCAN DESKTOP ICON DISPLAYS THE "IDLE" MESSAGE.



2. **If using the RF Power Measurement Kit** then refer to the RF Power Measurement Kit laminated card set.
 - a. Look in the upper, right corner of each card and find the card labeled **72** (1.5T Body Output).



The body RF output connection is no longer to the non-existent EFB unit, as the older RF Power Measurement Kit cards show, but instead to the J4 output on the rear of the SRFD2 Module. An RF adapter is provided in the RF Power Measurement Kit to connect between the HN J4 body RF output and the RF Power Measurement Kit 40 dB N-connector coupler.

- b. Configure the system as shown in the illustration on the card.
 - c. Confirm that the rotary attenuator is set to the correct position indicated on the card.
3. **If using the wattmeter or scope** (NOT the RF Power Measurement Kit) to measure power then refer to Appendix A (Alternate Equipment Setup) for the proper system body configuration.

3-1 Body Gain Pot Adjustment (continued)

4. Prepare the system to scan in Body mode per Table 3-1 or refer to Appendix A (Non-Proprietary protocol).

Note

If the system is operating with Release 9.X, CNV4 software, or its equivalent then it is no longer necessary to modify the **ia_rf1** and **ia_rf2** CV values. These will be set automatically when the **calmode** CV value is set equal to 5.

TABLE 3-1
SCAN PROTOCOL: BODY MODE

Note: This is the alternate proprietary procedure available for GE use and for sites with a valid Advanced Service Package Limited License.
Refer to Appendix A for the non-proprietary protocol.

A. **[New Pt]**

Id: **geservice<ENTER>**

Name: **rfi cals**

Weight (Lb.): **300<ENTER>**

Set Patient Protocols to **Service**.

B. At front enclosure:

Landmark in the Head area—remove any coils.

press **LANDMARK**.

press **MOVE TO SCAN**.

C. In the Patient Position Protocol field:

type **o.41.1<ENTER>** (o=Other, 41.1 =series) to load the body protocol

OR select **other** and select protocol **41** and select series **1**.

D. **[Save Series]**.

E. **[Research Operations]**.

[Setup Params]. Set TG to **50 [Done]**.

F. **[Research Operations]**.

[Display CVs]. Highlight CV Name and enter the following:

CV Name: **calmode<ENTER>**, CV Value: **5<ENTER>** (Dual Logamp Waveform).

Note

Skip the next two steps if the system has Release 9.X, CNV4, or equivalent software. The software will set the **ia_rf1** and **ia_rf2** CV values automatically once the **calmode** CV Value is set equal to 5.

CV Name: **ia_rf1<ENTER>**, CV Value: **32766<ENTER>** (sets 90° pulse full scale).

CV Name: **ia_rf2<ENTER>**, CV Value: **0<ENTER>** (turns off 180° pulse).

[Accept]

G. **[Research Operations]**

H. **[Download]**

5. Verify Body LED is illuminated on front of SRFD2 Module.
6. **[Manual Prescan] [Scan TR]**. Slowly increase TG to 200.

3-1 Body Gain Pot Adjustment (Continued)



Check the system center frequency NOW, especially if this is a new installation, and confirm that it consists of 8 digits and that it is correct. The RF amplifier will output a low and distorted RF waveform if center frequency is wrong due to a missing digit.

Note


If the RF Amplifier trips, adjust the Body Gain pot counter clockwise to reduce the RF output power and then increase the TG. Repeat this process until TG = 200 and the body RF output power is 16kW (72dBm).

Note

The Power Amplifier should be capable of 0.3dB to 0.4dB more than the required 72dBm RF Output. Improperly adjusted systems will trip if the gain is set above 72dBm creating a "Peak Power Fault".

7. Verify the Unblank LED on the front of the SRFD2 Module pulses on/off.

3-1 Body Gain Pot Adjustment (Continued)

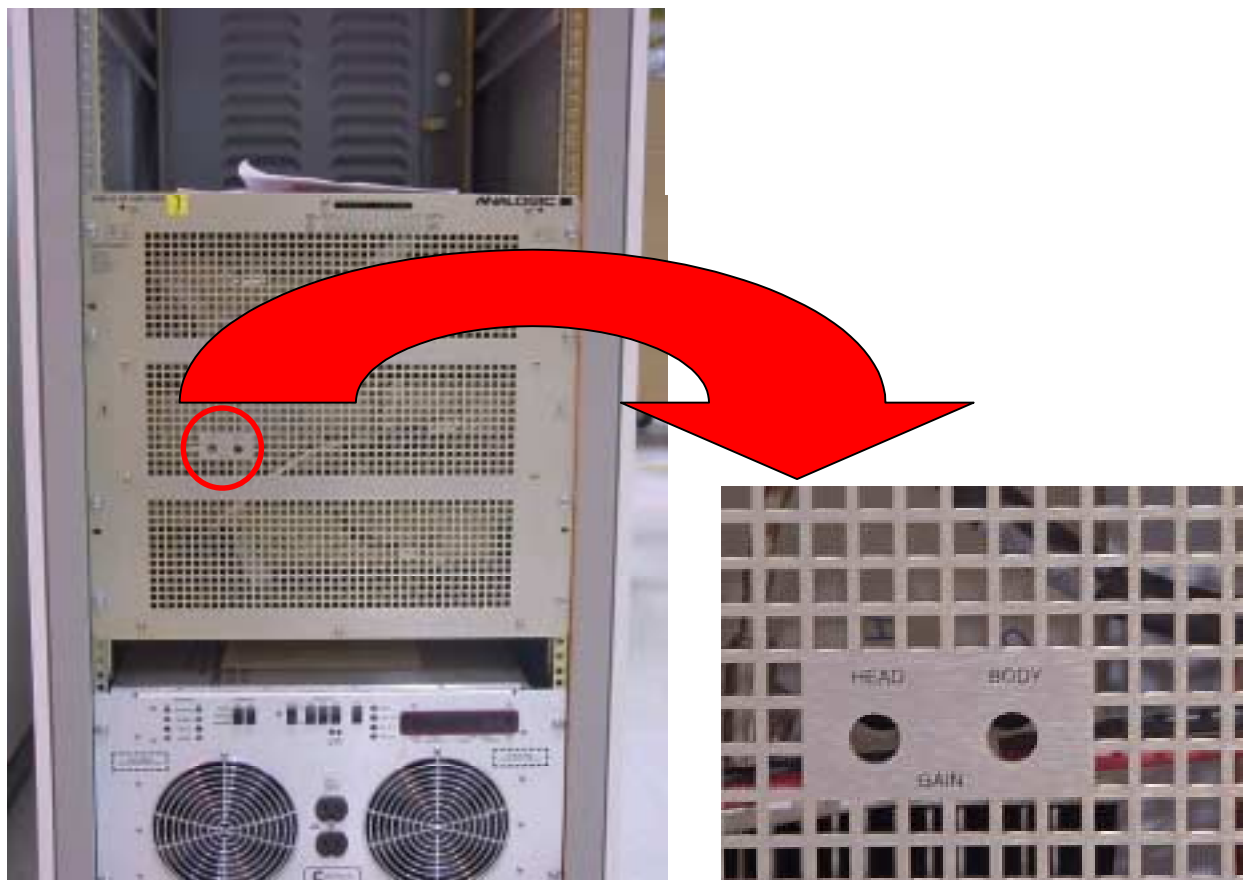

POT LOCATIONS ON SRFD2 MODULE
ILLUSTRATION 3-1

8. **If using the RF Power Measurement Kit**, refer to the specified number of divisions printed on the RF Power Measurement Kit card **72** (1.5T Body RF Output) needed to achieve 16kW of RF power output. Adjust the Body Gain pot on the front of the SRFD2 so that the RF waveform displayed on the scope meets, but does not exceed, this number of divisions. See Illustration 3-1.
9. **If using the wattmeter procedure** then read the wattmeter display and calculate the RF output power from the formula below. The dummy load and cable loss factor was determined from the procedure in Appendix F. Adjust the Body Gain pot on the front of the SRFD2 Module until the RF output calculated from the formula meets, and does not exceed, the 16kW (72dBm) specification. See Illustration 3-1.

RF Power Measurement (in watts) Using Wattmeter And Formula:
Wattmeter reading (in watts) X dummy load and cable loss factor

3-1 Body Gain Pot Adjustment (Continued)

10. If using the oscilloscope procedure (NOT the RF Power Measurement Kit) then read the peak voltage (V_{peak}) from the scope display and use the formula below or the [Power Calculator](#) Tool (located at E:\rf\power\pwrcalc.htm on the Service Methods CD-ROM) to calculate the RF output power. The dummy load and cable loss factors were determined from the procedure in Appendix F. The scope correction factor was determined in Appendix D. Adjust the Body Gain pot on the front of the SRFD2 until the RF output calculated from the formula meets, and does not exceed, the 16kW (72dBm) specification. See Illustration 3-1.

RF Power Measurement (in watts) Using Oscilloscope And Formula:
$\frac{\left(\frac{V_{\text{peak}}}{\text{scope correction factor}} \right)^2}{100} \times \text{dummy load and cable loss factor}$

11. Decrease TG to 0 (zero). **[Done]**.
12. Replace the RF Body heliax cable.

3-2 Head Gain Pot Adjustment



PROPERTY DAMAGE! PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, BY REMOVING ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

1. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the two **DANGER** messages on this page.



PERSONAL INJURY! PREVENT POSSIBLE RF BURNS WHEN DISCONNECTING HELIAX CABLES FROM J3 OR J4 ON THE SRFD2 BY VERIFYING THAT THE SYSTEM IS NOT MANUALLY PRESCANNING OR SCANNING. VERIFY THAT THE SCAN DESKTOP ICON DISPLAYS THE "IDLE" MESSAGE.



2. **If using the RF Power Measurement Kit** then refer to the RF Power Measurement Kit laminated card set.
 - a. Look in the upper, right corner of each card and find the card labeled **63** (1.5T Head Output).



The head RF output connection is no longer to the non-existent EFB unit, as the reference cards in some of the older kits show, but instead to the J3 output on the rear of the SRFD2.

- b. Configure the system as shown in the illustration on the card.
 - c. Confirm that the rotary attenuator is set to the correct position indicated on the card.
3. **If using the wattmeter or scope** (NOT the RF Power Measurement Kit) to measure power, refer to Appendix A (Alternate Equipment Setup) for the proper system head configuration.

3-2 Head Gain Pot Adjustment (Continued)

4. Prepare the system to scan in Head mode per Table 3-2 or refer to Appendix A (Non-Proprietary protocol).

Note

If the system is operating with Release 9.X, CNV4 software, or its equivalent then it is no longer necessary to modify the **ia_rf1** and **ia_rf2** CV values. These will be set automatically when the **calmode** CV value is set equal to 5.

TABLE 3-2
SCAN PROTOCOL: HEAD MODE

Note: This is the alternate proprietary procedure available for GE use and for sites with a valid Advanced Service Package Limited License.
Refer to Appendix A for the non-proprietary protocol.

- A. **[New Series]**
At Patient Protocols – select **other**.
- B. In the protocol field, type **o.41.2<ENTER>** (o=Other, 41.2 =series) to load the head protocol
OR select **[o.41] [Series 2] [Accept]**.
- C. **[OK]** (if required).
- D. **[Save Series]**
- E. **[Prepare to Scan]**
- F. **[Research Operations]**.
[Setup Params]. Set TG to **50 [Done]**.
- G. **[Research Operations]**.
[Display CVs]. Highlight CV Name and enter the following:
CV name: **calmode<ENTER>, 5<ENTER>** (Dual Logamp Waveform).

Note

Skip the next two steps if the system has Release 9.X, CNV4, or equivalent software. The software will set the **ia_rf1** and **ia_rf2** CV values automatically once the **calmode** CV Value is set equal to 5.

CV name: **ia_rf1<ENTER>, 32766<ENTER>** (sets 90° pulse full scale).
CV name: **ia_rf2<ENTER>, 0<ENTER>** (turns off 180° pulse).
[Accept].

- H. **[Research Operations]**
- I. **[Download]**.


5. Verify Head LED is illuminated on front of SRFD2 Module.



The Head Gain pot is directly affected by changes to the Body Gain pot. Adjust the Head Gain pot ONLY after completing final adjustment of the Body Gain pot.

6. **[Manual Prescan] [Scan TR]**. Increase TG to 200.
7. Verify the Unblank LED on the front of the SRFD2 Module pulses on/off.

3-2 Head Gain Pot Adjustment (Continued)

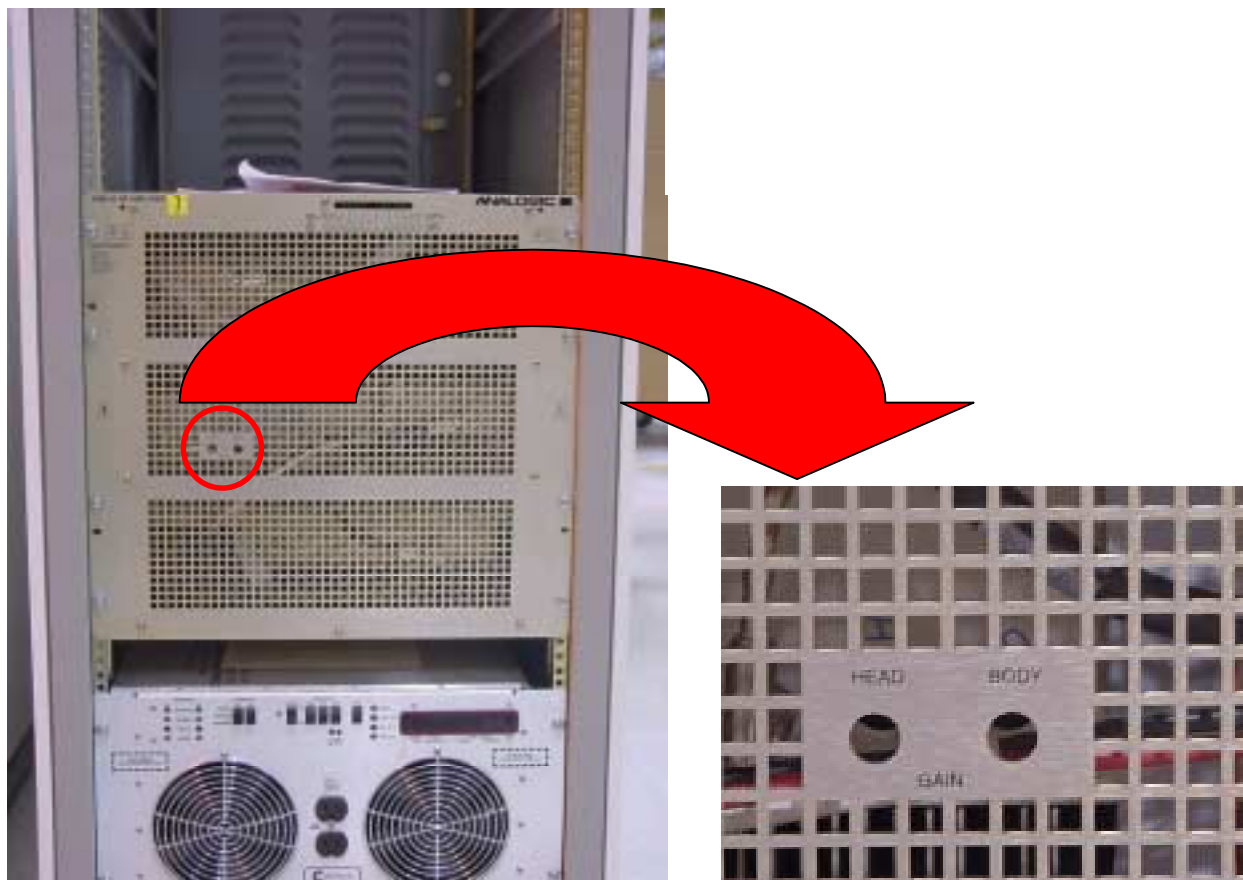

POT LOCATIONS ON RFI MODULE
ILLUSTRATION 3-4

8. **If using the RF Power Measurement Kit** then refer to the specified number of divisions printed on the RF Power Measurement Kit card **63** (1.5T Head RF Output) needed to achieve 2kW of RF power output. Adjust the Head Gain pot on the front of the SRFD2 so that the RF waveform displayed on the scope meets, but does not exceed, this number of divisions. See Illustration 3-4.
9. **If using the wattmeter procedure** then read the wattmeter display and calculate the RF output power from the formula below. The cable loss factor was determined from the procedure in Appendix F. Adjust the Head Gain pot on the front of the SRFD2 until the RF output calculated from the formula meets, and does not exceed, the 2kW (63dBm) specification. See Illustration 3-4.

RF Power Measurement (in watts) Using Wattmeter And Formula:
Wattmeter reading (in watts) X cable loss factor


10. If using the oscilloscope procedure (NOT the RF Power Measurement Kit) then read the peak voltage (V_{peak}) from the scope display and use the formula below or the [Power Calculator](#) Tool (located at E:\rf\power\pwrcalc.htm on the Service Methods CD-ROM) to calculate the RF output power. The dummy load and cable loss factor was determined from the procedure in Appendix F. The scope correction factor was determined in Appendix D. Adjust the Head Gain pot on the front of the SRFD2 module until the RF output calculated from the formula meets, and does not exceed, the 2kW (63dBm) specification. See Illustration 3-4.

RF Power Measurement (in watts) Using Oscilloscope And Formula:
$\frac{\left(\frac{V_{\text{peak}}}{\text{scope correction factor}} \right)^2}{100} \times \text{dummy load and cable loss factor}$

11. Decrease TG to 0 (zero). **[Done] [End Exam]**.
12. Verify scan desktop icon is displaying the “Idle” message.
13. Replace the RF Head Heliax cable.
14. Proceed to **Section 4 - System Restoration**.

4 - SYSTEM RESTORATION

1. Verify the system is not scanning and the scan desktop icon is displaying the "Idle" message.
2. See Illustration 4-1. On the front panel of the SSM place the:
 - 2 (two) power MONITOR switches (A and B) to the bottom Normal position.
 - 3 (three) DRV FLT switches to the top EN (enable faults) position.

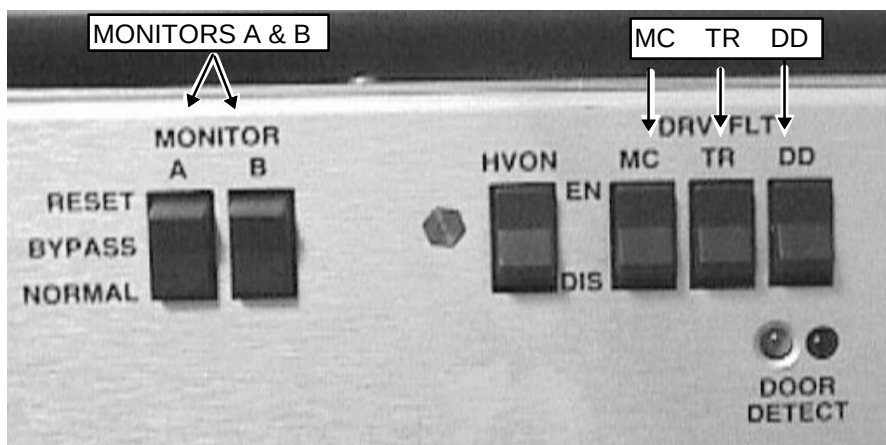

SSM FRONT PANEL SWITCHES ENABLED
ILLUSTRATION 4-1


3. Remove all test equipment.
4. Verify the Body Heliax cable is connected.
5. Verify the Head Heliax cable is connected.
6. Re-install the cabinet covers.
7. Successfully complete one body scan.
8. Successfully complete one head scan.

APPENDIX A — ALTERNATE EQUIPMENT SETUPS & NON-SERVICE PROTOCOLS

Use this section if the RF Power Measurement Kit is NOT going to be used. This section contains the setup information and protocols necessary to measure the body and head RF output power using either a wattmeter or oscilloscope. If the wattmeter is not used it is strongly recommended that an oscilloscope with a bandwidth of 300 MHz or greater be used for making the RF output power measurements. A 100 MHz oscilloscope can be used to make these measurements **provided** that the oscilloscope input channel to be used has been properly characterized and the resulting correction factor is known. The oscilloscope must be characterized and the correction factor must be known to account for measurement error at 63.86 MHz error due to bandwidth limitations of the 100 MHz oscilloscope.

TABLE A-1
EQUIPMENT REQUIRED IF NOT USING THE RF POWER MEASUREMENT KIT

Item	Description	Part Number
1.	RF Test Cable Kit	46-255816G1
2.	100 MHz Scope (equivalent or greater)	46-183029P61
3.	Digital Multimeter	46-194427P49
4.	50 ohm, 200 Watt, 30 dB Attenuator (dummy load)	46-317724P14
5.	Screwdriver, slotted (long & narrow)	—
6.	Wattmeter Kit (optional)	46# not supplied
7.	DMM Test Cable – BNC to Banana (optional)	2239139

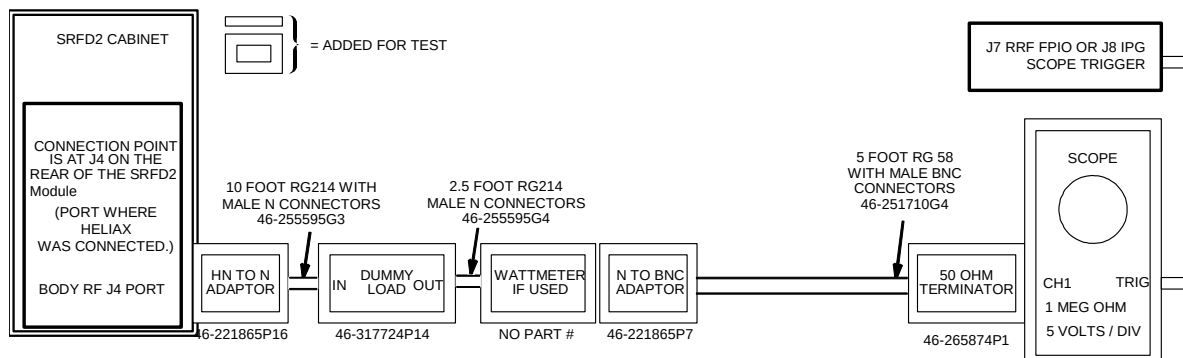

ALTERNATE BODY EQUIPMENT SETUP FOR WATTMETER OR SCOPE MEASUREMENT
ILLUSTRATION A-1

Note

If the wattmeter is not available then bypass it's connection in the circuit using an N-female to N-female connector 46-265875P2.

Note

If the system is operating with Release 9.X, CNV4 software, or its equivalent then it is no longer necessary to modify the ia_rf1 and ia_rf2 CV values. These will be set automatically when the calmode CV value is set equal to 5.

NON-SERVICE BODY PROTOCOL

PATIENT REGISTER .
PATIENT INFORMATION

Patient Id **geservice**
Patient Name **body rf**
Weight (Lb) **300 ~ IMPORTANT**
Landmark **[Landmark]**
[>] [Sternal Notch]

PATIENT PROTOCOLS **[Patient Position]**

PATIENT POSITION

Patient Position **[>] [Supine]**
Patient Entry **[>] [Head First]**
Coil **[...] [Body] [Accept]**

IMAGING PARAMETERS

Plane **[>] [Axial]**
Mode **[>] [2D]**
Pulse Seq **[...] [Spin Echo]**
[Accept]
Imaging Options none (default)
Psd Name **cal**
Protocol no entry

SCAN TIMING

* of Echoes 1 (default)
TE **[25]**
TR **[55]**

SCANNING RANGE

FOV **[24]**
Slice Thickness **[5]**
Spacing **0**
Start **0**
End **0**
Slices 1 (default)
L/R Center 0 (default)
P/A Center 0 (default)
Table Delta 0.00 (default)

ACQUISITION TIMING

Freq **[256]**
Phase **[128]**
NEX **[2]**
Freq Dir **[>] [A/P]**
Auto Center Freq **[>] [Peak]**
(lowest window) **[Save Series]**

[Research Operations] [Display CVs]

Modify the following: (if 9.X, CNV4, or equivalent only modify **calmode**)


calmode **5**
ia_rf1 **32766** (90° maximum)
ia_rf2 **0** (180° minimum)
[Accept]

[Research Operations] [Setup Params]

Set TG **50 [Done]**

[Research Operations] [Download]

[Prepare to Scan]

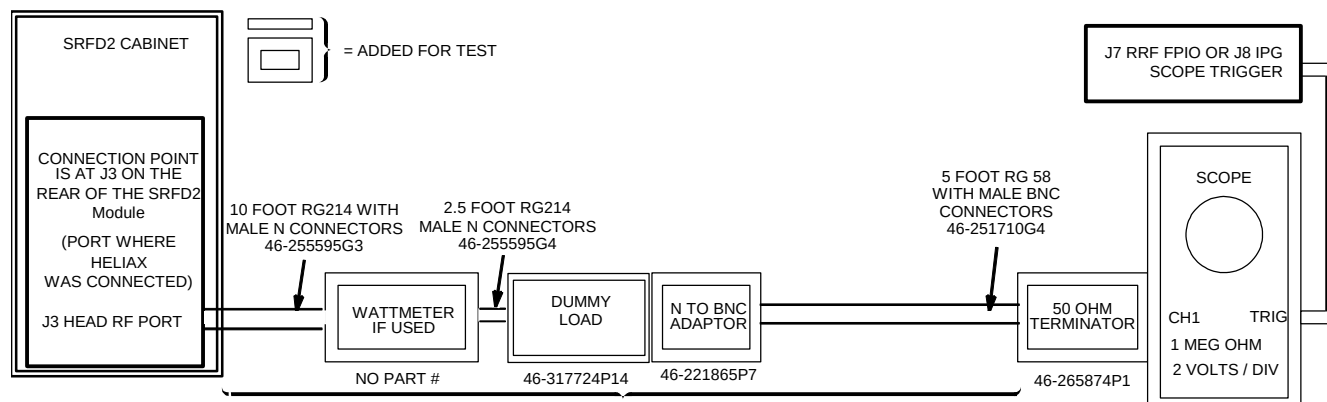

ALTERNATE HEAD EQUIPMENT SETUP
ILLUSTRATION A-2

Note

If the wattmeter is not available then bypass it's connection in the circuit using an N-female to N-female connector 46-265875P2.

Note

If the system is operating with Release 9.X, CNV4 software, or its equivalent then it is no longer necessary to modify the ia_rf1 and ia_rf2 CV values. These will be set automatically when the calmode CV value is set equal to 5.

NON- SERVICE HEAD PROTOCOL

<u>PATIENT REGISTER .</u>	[New Pt]	<u>SCANNING RANGE</u>	
<u>PATIENT INFORMATION</u>		FOV	[24]
Patient Id	geservice	Slice Thickness	[5]
Patient Name	head rf	Spacing	0
Weight (Lb)	300 - IMPORTANT	Start	0
	[Landmark]	End	0
Landmark	[>] [Sternal Notch]	# Slices	1 (default)
		L/R Center	0 (default)
<u>PATIENT PROTOCOLS</u>	[Patient Position]	P/A Center	0 (default)
		Table Delta	0.00 (default)
<u>PATIENT POSITION</u>		<u>ACQUISITION TIMING</u>	
Patient Position	[>] [Supine]	Freq	[256]
Patient Entry	[>] [Head First]	Phase	[128]
Coil	[...] [Head] [Accept]	NEX	[2]
<u>IMAGING PARAMETERS</u>		Freq Dir	[>] [A/P]
Plane	[>] [Axial]	Auto Center Freq	[>] [Peak]
Mode	[>] [2D]	(lowest window)	[Save Series]
Pulse Seq	[...] [Spin Echo]	[Research Operations] [Display CVs]	
	[Accept]	Modify the following: (if 9.X, CNV4, or equivalent only modify calmode)	
Imaging Options	none (default)	calmode	5
Psd Name	cal	ia_rf1	32766 (90° maximum)
Protocol	no entry	ia_rf2	0 (180° minimum)
<u>SCAN TIMING</u>			[Accept]
* of Echoes	1 (default)	[Research Operations] [Setup Params]	
TE	[25]	Set TG	50 [Done]
TR	[55]	[Research Operations] [Download]	
		[Prepare to Scan]	

APPENDIX B — CALIBRATION REQUIREMENTS AFTER HARDWARE REPLACEMENT

TABLE B-1
RE-CALIBRATION REQUIREMENTS AFTER HARDWARE REPLACEMENT

Hardware Replaced	Required Calibration
1. Exciter Board ((U)CERD)	1. Entire procedure
2. SRFD2 Module	2. Entire procedure
3. RRF Master Exciter	3. Entire procedure
4. RRF Multiplexer	4. Entire procedure

APPENDIX C — SOFTWARE AND FORMULAE FOR MISCELLANEOUS RF CONVERSIONS

Software are described and formulae are presented in this section that will allow the FE to perform various conversions and calculations.

C-1 RF Conversion Software

C-1-1 RF Power Calculator

The first tool, the RF Power Calculator, is a UNIX-based application available on the LX system under the **Toolbelt icon**, **[Utilities]**, **[RF Calculator]**, **[Start]**. Watts to dBm, dBm to Watts, relative volts to relative power, and relative power to relative volts conversions can be done with the RF Power Calculator tool.

C-1-2 Power Calculator Tool

The second tool, the [Power Calculator](#), is an HTML document that is available from the 8.X Service Methods CD-ROM. This tool will run from either the laptop PC or MR host computer. The tool can manually be started from the laptop PC by selecting **Start**, **Run**, and then entering the following file pathname in the “Open” entry box: E:\rf\power\pwrcalc.htm. Next, select **OK**. The tool is divided into 3 sections. The first section allows one to calculate RF power from measured peak voltage when using a 100 MHz oscilloscope and a properly characterized RF cables kit and dummy load. The second section provides Vpp to dBm, dBm to Watts, and dBm to Vpp conversion capability. The third section provides short instructions and a means for comparing the attenuation in dB reported on each RF Power Measurement Kit Card with the measured Magnitude Squared Attenuation Factor reported by the Attenuation Test tool that is described in **APPENDIX D – DUMMY LOAD AND CABLES CALIBRATION**.

C-2 RF Conversion Formulae

These are provided to assist the experienced FE with RF power measurement and troubleshooting. The RF conversions below are presented in two different formats. The first format is the actual mathematical formula. The second format assumes that you are entering the data into the LX system calculator exactly as you read it (Access the LX system calculator by right mouse clicking in the background and then selecting the calculator from the Root menu.).

VP-P to dBm Calculation:

$$\text{dBm} = 20 \log \left(\frac{V_{pp}}{0.632} \right)$$

VP-P **[+]** 0.632 **[=]** **[LOG]** **[*]** 20 **[=]** dBm.

APPENDIX C — SOFTWARE AND FORMULAE FOR MISCELLANEOUS RF CONVERSIONS (CONTINUED)

dBm to VP-P Calculation:

$$V_{pp} = 10^{\left(\frac{dBm}{20}\right)} \times 0.632 \quad \text{or} \quad V_{pp} = \text{antilog} \left(\frac{dBm}{20} \right) \times 0.632$$

dBm [+] 20 [=] [INV LOG] [*] 0.632 [=] VP-P.

Example: Using 3.65 dBm, the Scope Power reading listed on a particular Kit card, and the LX system calculator to determine the equivalent Vpp:

3.65 [+] 20 [=] (0.1825) [INV LOG] (1.5222991) [*] 0.632 [=] (0.9620931) VP-P.

dBm to Watts Calculation:

$$\text{Watts} = 10^{\left(\frac{dBm}{10}\right)} \times 0.001 \quad \text{or} \quad \text{Watts} = \text{antilog} \left(\frac{dBm}{10} \right) \times 0.001$$

dBm [+] 10 [=] [INV LOG] [=] [*] 0.001 [=] Total Watts.

C-3 1.5T Division Conversions

The 1.5T division conversion below may be useful for the FE using the RF Power Measurement Kit and attempting to measure a signal that requires more resolution than the "division" unit provides. These examples are provided as a general reference and are not intended to be used to perform the power monitor check.

1.5T Division Conversion example:

- 1.5T sites using the 1.5T Scope Calibrator and 100 MHz Oscilloscope equipment:

$$\frac{1.00 \text{ VPP (4.0 dBm)}}{\text{Set to 8 Major Div}} = \frac{0.125 \text{ VPP each Major Division}}{5 \text{ Minor Divisions}} = 0.025 \text{ VPP each Minor Division}$$

$$0 \text{ dBm} = \frac{0.632 \text{ VPP}}{0.025 \text{ VPP each Minor Division}} = 25.28 \text{ Minor Divisions}$$

25.28 Minor Divisions = 5 Major Divisions and .28 of a Minor Division

APPENDIX D — CHARACTERIZATION FOR SCOPES WITH < 300 MHZ BANDWIDTH

WARNING!

THIS PROCEDURE REQUIRES THAT THE SCOPE IN USE HAS A BANDWIDTH OF AT LEAST 100 MHZ. DO NOT ATTEMPT TO PERFORM RF POWER MEASUREMENTS WITH AN OSCILLOSCOPE THAT HAS A BANDWIDTH LESS THAN 100 MHZ. SIGNIFICANT MEASUREMENT ERRORS CAN RESULT.

The scope correction factor for the scope channel in use must be known for oscilloscopes with a bandwidth < 300 MHz. One may assume that the correction factor is 1 (unity) for oscilloscopes with a bandwidth ≥ 300 MHz. Do not assume a correction factor for an input channel and do not assume that all the scope input channels have the same correction factor. Each channel must be characterized separately. There are three ways that this can be done.

D-1 Scope Channel Characterization Done During Calibration

One way that the scope channel characterization can be accomplished is to request this information from the metrology lab the next time the scope is sent back for calibration. Attach a tag in a conspicuous location on the scope before it is sent back for calibration requesting this information. Ask that the percent of signal loss at 63.86 MHz (1.5T) and 42 MHz (1.0T) be measured at the 5 volt/div and 2 volt/div settings (for body and head measurements, respectively) for both scope channels and that this information be included with the returned scope. The metrology lab can usually do this upon request and without charge during the calibration process.

D-2 Scope Channel Characterization Done By Comparison

A 400 MHz scope exhibits negligible measurement error at 63.86 MHz or 42 MHz. If a 400 MHz scope is available then it can be compared side-by-side on site with a scope of < 300 MHz bandwidth. Be sure the scope inputs have a 1Mohm impedance and then are terminated with a 50 ohm feed-through adapter before making the measurements. Be sure to characterize all the inputs that are to be used. The signal loss at 63.86 MHz (1.5T) and 42 MHz (1.0T) for each channel at the 5 volt/div and 2 volt/div settings can be determined from the ratio of signal measured on the 100 MHz scope to that seen on the 400 MHz scope

$$\frac{100 \text{ MHz scope measured Vpp signal}}{400 \text{ MHz scope measured Vpp signal}} = \text{scope correction factor}$$


. As long as the same cables are used the loss can also be determined for each cable.

D-3 Scope Characterization Using the RF Power Measurement Kit Signal Reference

The RF Power Measurement Kit contains a 63.86 MHz, 4dBm (1 Vpp) calibrated signal reference that can be used to characterize the scope input channels. Be aware that these values are accurate only for 1.5T systems.

APPENDIX D (CONTINUED)

1. Connect the calibrated reference to the scope as in Illustration D3-1 below.

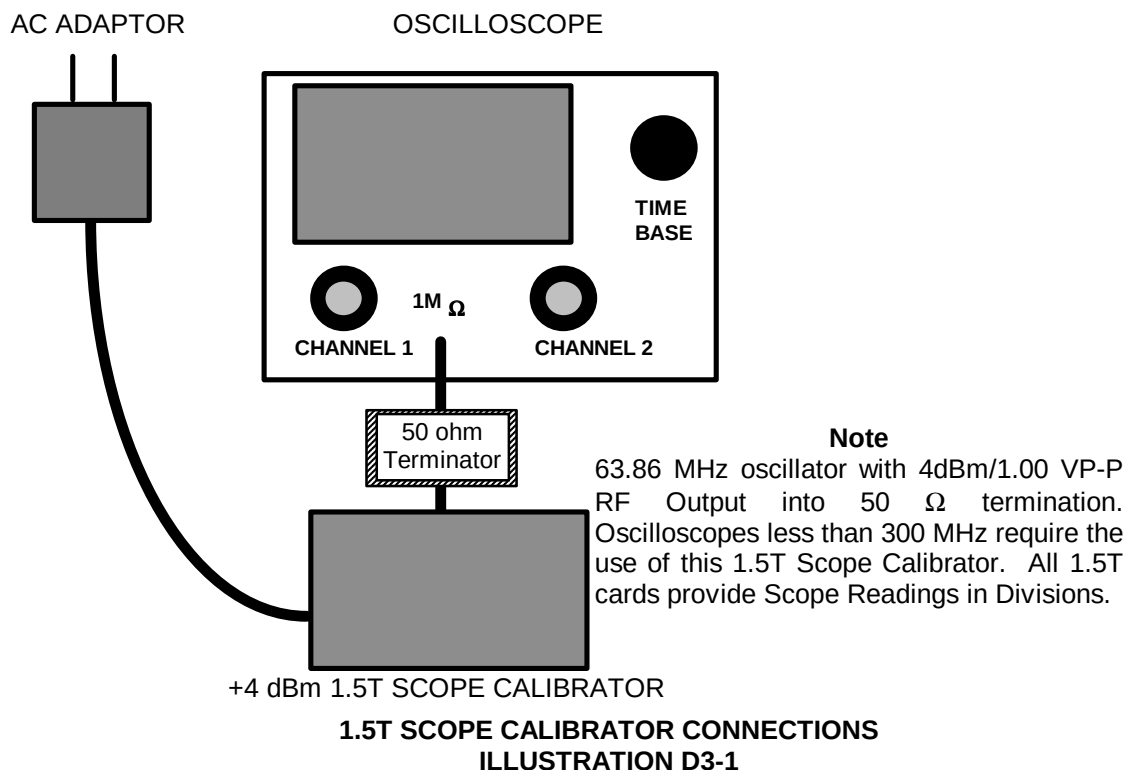

2. Set the scope amplitude for the channel to be checked to 0.2 volts/div.
3. Set the sweep rate to 5 nsec.
4. Confirm that the channel input impedance is set for 1Mohm and not 50 ohms.
5. Confirm that the scope bandwidth limit switch is not enabled.
6. Confirm that the "Uncal" knob is not rotated and the channel is not uncalibrated.
7. Carefully measure the peak to peak voltage from the calibrated source. The measured value should be 1 Vpp or less.
8. Calculate the correction factor (it should be less than 1):

$$\frac{V_{pp} \text{ measured}}{1 V_{pp}} = \text{correction factor}$$

9. Repeat this process and derive the correction factor for the other scope input channel(s).
10. Note these correction factors for each channel on a piece of tape affixed to the scope for later reference.

APPENDIX E — (U)CERD EXCITER OR RRF MASTER EXCITER RF SIGNAL OUTPUT

A minimum RF signal level from the (U)CERD Exciter or RRF (Remote RF) Master Exciter is required or else it will be impossible to meet the specified 16kW body and 2kW head RF power output levels. Use this section to verify that the (U)CERD Exciter or RRF Master Exciter is outputting at least the minimum amount of power required to the SRFD2. See Table E-1 for needed tools and test equipment.

TABLE E-1
ITEMS NEEDED TO CHECK (U)CERD RF SIGNAL OUTPUT

Item	Description	Part Number
1.	RF Test Cable Kit	46-255816G1
2.	Cannon to BNC-male Test Cable (included in item 4 below)	46-301549P6
3.	SMB-female to BNC-male Test Cable (included in item 4 below)	46-301549P5
4.	TPS RF Connector/Adapter and Cable Test Kit (optional)	46-301927G1
5.	RF Power Measurement Kit (recommended)	46-317724G1 or G2
6.	100 MHz Scope (equivalent or greater)	46-183029P61


Note

Please read the following if you do not have an RF Power Measurement Kit and you are going to attempt to measure the exciter output directly with a 100 MHz oscilloscope. Be aware that, due to oscilloscope bandwidth limitations, an exciter RF voltage reading taken directly from a 50 ohm terminated 100 MHz oscilloscope input channel can be as much as 16% less than the actual exciter RF output voltage level. This error is accounted for and is not a problem when using the oscilloscope in conjunction with the RF Power Measurement Kit. Reading *reasonably* accurate 63.86 MHz voltage levels directly from a 100 MHz oscilloscope *without* using the RF Power Measurement Kit requires that the correction factor of the scope channel being used is known. This is so that the approximate true RF voltage level can be calculated. The FE can request that the 63.86 MHz correction factor for each channel be determined and reported by the metrology lab the next time the oscilloscope is returned for calibration or, if the FE has access to the 4dBm (1 Vpp) oscilloscope calibrator tool used in the RF Power Measurement Kit, he can determine this correction factor himself. See Appendix D for more information concerning the derivation of the correction factor. Divide the Vpp value read from the oscilloscope by the correction factor to determine the approximate true RF voltage level value. The simple formula is as follows:

$$\frac{\text{Measured V}_{\text{peak to peak}}}{\text{Correction factor}} = \text{Approximate True RF Voltage Value}$$

E-1 Exciter or RRF Master Exciter RF Output Measurement Using an Oscilloscope

1. Connect the oscilloscope to the system hardware as per Illustration E1 -1. Confirm that the oscilloscope is properly configured, that the channel 1 vertical Volts/div variable control is fully CCW, and that the bandwidth limit button on the oscilloscope is not selected.

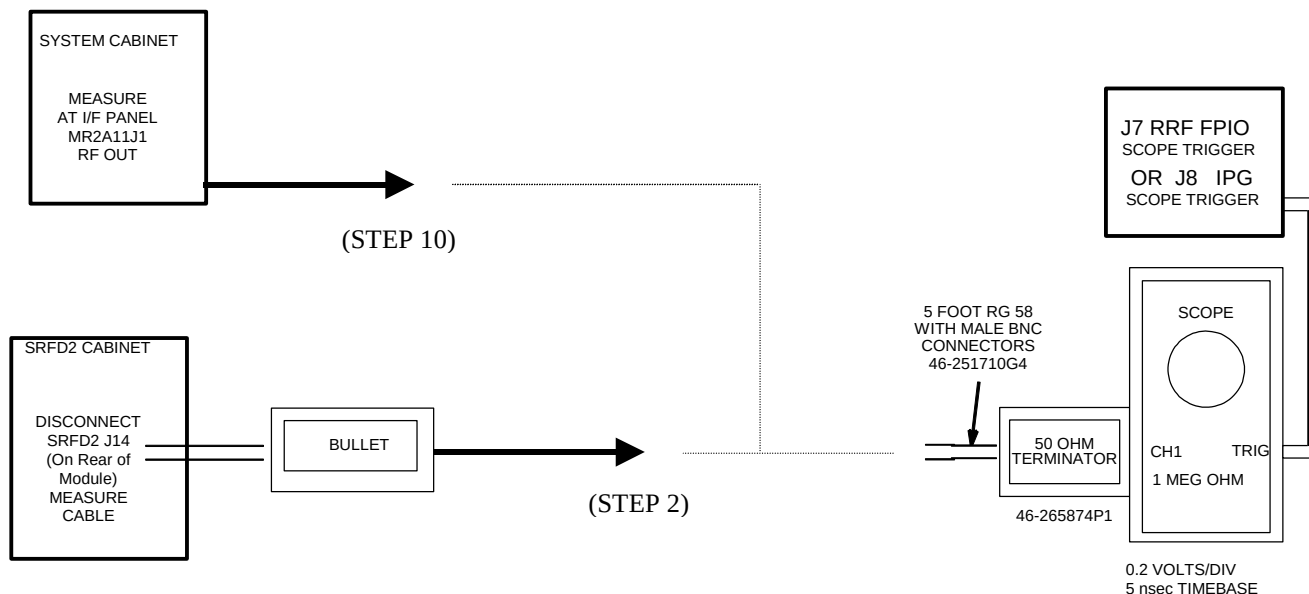

RF SIGNAL MEASUREMENTS USING A ≥ 300 MHZ OR CHARACTERIZED 100 MHZ SCOPE ILLUSTRATION E1-1

2. Ensure that channel 1 is terminated with a 50 ohm feed-through terminator. Disconnect the cable going to J14 on the rear of the SRFD2 module and connect it to the input of the 50 ohm terminator connected to channel 1 on the scope. A short (less than 5 ft.) length of coax and a BNC bullet (female BNC/female BNC) adapter can be connected to the cable removed from J14 to route the signal to the scope.
3. Prepare the system to scan using the proprietary protocol listed in Table 2-3.
4. **[Manual Prescan] [Scan TR]**. Increase TG to 200.
5. If using a scope with a bandwidth of 300 MHz or greater then read the peak to peak voltage from the scope face and confirm that it is ≥ 0.632 Vpp (0 dBm).
6. If using a < 300 MHz scope read the peak to peak voltage from the scope face and then, *using the correction factor for the scope channel in use*, calculate the approximate true RF voltage value by using the formula below. Confirm that the approximate true RF voltage value is ≥ 0.632 Vpp (0 dBm).

$$\frac{\text{Measured } V_{\text{peak to peak}}}{\text{Correction factor}} = \text{Approximate True RF Voltage Value}$$

7. Decrease TG to 0 (zero). **[Done]**.

E-1 Exciter RF Output Measurement Using an Oscilloscope (Continued)

8. Reconnect the cable removed from J14 on the rear of the SRFD2.
9. If the measured voltage meets the specification then the Exciter is working properly. The functionality of the SRFD2 should be checked again. No need to continue this procedure.
10. Disconnect the cable going to J1 on the rear of the System Cabinet and connect to J1 a known good coax test cable 5 ft. in length. Route the other end of the test cable to the input of the 50 ohm feed-through terminator connected to channel 1 on the scope.
11. **[Manual Prescan] [Scan TR]**. Increase TG to 200.
12. If using a 300 MHz bandwidth or greater scope then read the peak to peak voltage from the scope face and confirm that it is 0.893 to 1.12 Vpp (3.0 dBm to 5.0 dBm). A voltage level exceeding the 1.12 Vpp (5.0 dBm) typical upper limit is generally not a cause for concern.
13. If using a 100 MHz bandwidth scope read the peak to peak voltage from the scope face and then, *using the correction factor for the scope channel in use*, calculate the approximate true RF voltage value by using the formula in step 6 and confirm that the approximate true RF voltage value is ≥ 0.893 Vpp (3.0 dBm).
14. Decrease TG to 0 (zero). **[Done]**.
15. Reconnect the cable removed from J1 on the rear of the System Cabinet.
16. Check for poor connections in the cabling or connectors between the System Cabinet and RF Cabinet if the voltage measured at J14 is too low but the voltage at J1 on the rear of the System Cabinet meets specification.
17. **If the system has a (U)CERD:** Remove connector from J109 on the bottom, front of the (U)CERD.
 - a. Insert the cannon-to-BNC test cable from the "TPS RF Connector/Adapter and Cable Test Kit" into the top-most female plug inside the J109 connector on the (U)CERD, marked EXC RF OUT.
 - b. Connect the BNC end of the test cable to the 50 ohm feed-through connector on scope channel 1.
18. **If the system has an RRF:** Locate the Master Exciter Output J8 SMB port on the Multiplexer Board located on the rear of the RRF.
 - a. Remove the cable connected to the Master Exciter Output J8 connector and set it aside.
 - b. Insert the SMB-female to BNC-male test cable from the "TPS RF Connector/Adapter and Cable Test Kit" into the J8 Master Exciter Output connector.

E-1 Exciter or RRF Master Exciter RF Output Measurement Using an Oscilloscope (Continued)

- c. Connect the BNC end of the test cable to the 50 ohm feed-through connector on scope channel 1.
19. **[Manual Prescan] [Scan TR].** Increase TG to 200.
20. If using a 300 MHz bandwidth or greater scope then read the peak to peak voltage from the scope face and confirm that it is 0.893 to 1.12 Vpp (3.0 dBm to 5.0 dBm). A voltage level exceeding the 1.12 Vpp (5.0 dBm) typical upper limit is generally not a cause for concern.
21. If using a 100 MHz bandwidth scope read the peak to peak voltage from the scope face and then, *using the correction factor for the scope channel in use*, calculate the approximate true RF voltage value by using the formula in step 6 and confirm that the approximate true RF voltage value is ≥ 0.893 Vpp (3.0 dBm).
22. Decrease TG to 0 (zero). **[Done].**
23. **If the system has a (U)CERD:** Remove the cannon-to-BNC test cable from the top-most plug of J109 on the (U)CERD and reconnect the cable that was previously removed from J109.
24. **If the system has an RRF:** Remove the SMB-female to BNC-male test cable from the Master Exciter Output J8 connector on the Multiplexer Board and reconnect the cable that was previously removed.
25. **If the system has an RRF:** Trace the cable connecting to the Master Exciter Input J10 on the Multiplexer Board over to the J6 RF Out on the Master Exciter Board.
 - a. Remove the cable at J6 RF Out on the Master Exciter Board and connect the SMB-female to BNC-male test cable to this port.
 - b. Repeat steps 19 – 22 above.
 - c. Remove the SMB-female to BNC-male test cable from the J6 RF Out connector and reconnect the cable that was previously removed.
26. Only consider replacing the Exciter or RRF Master Exciter Boards if the output is below the specification **AND** the body RF output power cannot be adjusted to meet the 16 kW specification. Otherwise, check the RF cabling, connectors, or any hardware in between.

E-2 Exciter or RRF Master Exciter RF Output Check Using The RF Power Measurement Kit


More accurate measurements are made by using the RF Power Measurement Kit and a 100 MHz oscilloscope. The kit will compensate for any amplitude error due to the bandwidth limitation associated with using a 100 MHz scope. The 1.5T Scope Calibrator provides a 63.86 MHz sine wave at 1.00 VPP (4dBm).

1. Verify that the Bandwidth Limit button is not selected on the oscilloscope.
2. Verify Channel 1 is set to 1 Mega ohm input termination.

Note

Use the 50 ohm terminator supplied with the RF Power Measurement Kit and set Channel 1 oscilloscope termination selection to 1Mega ohm.

3. Connect the 1.5T Scope Calibrator via the 50 ohm feed-through terminator to the oscilloscope input. See Illustration E2-1 or the laminated card included in the RF Power Measurement Kit that is labeled, in the upper right-hand corner, CAL.

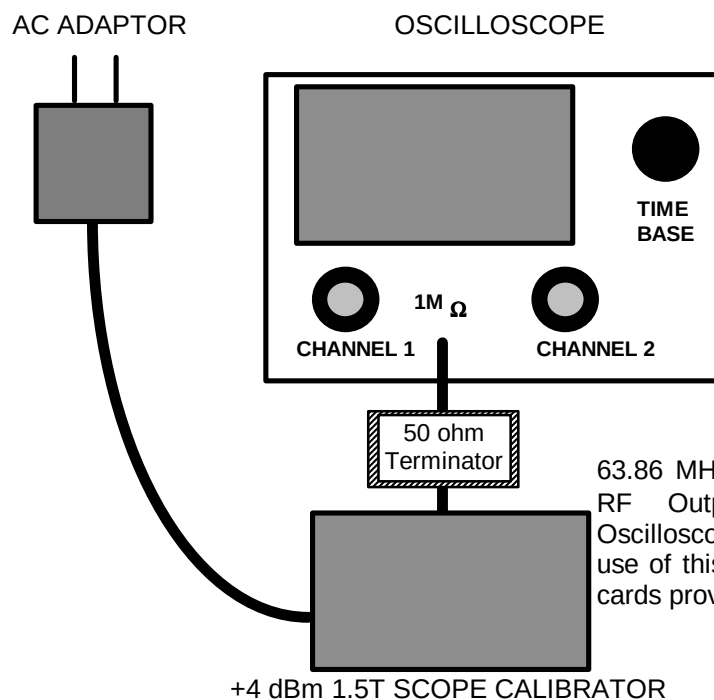

Note

63.86 MHz oscillator with 4dBm/1.00 VP-P RF Output into 50 Ω termination. Oscilloscopes less than 300 MHz require the use of this 1.5T Scope Calibrator. All 1.5T cards provide Scope Readings in Divisions.

1.5T SCOPE CALIBRATOR CONNECTIONS
ILLUSTRATION E2-1

4. Adjust the channel 1 vertical Volts/div control until amplitude of the waveform is slightly greater than 8 divisions.
5. Adjust the channel 1 vertical Volts/div variable control to achieve exactly 8 divisions.
6. Remove 1.5T Scope Calibrator. Do not remove the 50 ohm feed-through terminator attached to the oscilloscope channel 1.

E-2 Exciter RF Output Check Using The RF Power Measurement Kit (Continued)

7. Disconnect the cable going to J14 on the rear of the SRFD2 and connect it to the input of the 50 ohm terminator connected to channel 1 on the scope. A short (less than 5 ft.) length of coax and a BNC bullet (female BNC/female BNC) adapter can be connected to the cable removed from J14 to route the signal to the scope.
8. Prepare the system to scan using the proprietary protocol listed in Table 2-3.
9. **[Manual Prescan] [Scan TR].** Increase TG to 200.
10. Confirm that the output from the J14 cable is ≥ 25.28 **minor divisions**.

Note

Understand that each of the 8 divisions on the scope face is composed of 5 minor divisions. Since $8 \times 5 = 40$ total minor divisions, count the number of individual minor divisions that the RF waveform spans. Since the input signal is small, minor divisions are used in place of divisions so that a more accurate scope measurement can be made. This is discussed in detail in Appendix C.

11. Decrease TG to 0 (zero). **[Done].**
12. Reconnect the cable removed from J14 on the rear of the SRFD2.
13. If the measured voltage meets the specification then the Exciter is working properly. The functionality of the SRFD2 should be checked again. No need to continue this procedure.
14. Disconnect the cable going to J1 on the rear of the System Cabinet and connect to J1 a known good coax test cable 5 ft. in length. Route the other end of the test cable to the input of the 50 ohm feed-through terminator connected to channel 1 on the scope.
15. **[Manual Prescan] [Scan TR].** Increase TG to 200.
16. Read the peak to peak size of the waveform in divisions and then confirm that it is ≥ 35.7 **minor divisions**.
17. Decrease TG to 0 (zero). **[Done].**
18. Reconnect the cable removed from J1 on the rear of the System Cabinet.
19. Check for poor connections in the cabling or connectors between the System Cabinet and RF Cabinet if the voltage measured at J14 is too low but the voltage at J1 on the rear of the System Cabinet meets specification.
20. **If the system has a (U)CERD:** Remove connector from J109 on the bottom, front of the (U)CERD.

E-2 Exciter or RRF Master Exciter RF Output Check Using The RF Power Measurement Kit (Continued)

- a. Insert the cannon-to-BNC test cable from the “TPS RF Connector/Adapter and Cable Test Kit” into the top-most female plug inside the J109 connector on the (U)CERD, marked EXC RF OUT.
 - b. Connect the BNC end of the test cable to the 50 ohm feed-through connector on scope channel 1.
21. **If the system has an RRF:** Locate the Master Exciter Output J8 SMB port on the Multiplexer Board located on the rear of the RRF.
 - a. Remove the cable connected to the Master Exciter Output J8 connector and set it aside.
 - b. Insert the SMB-female to BNC-male test cable from the “TPS RF Connector/Adapter and Cable Test Kit” into the J8 Master Exciter Output connector.
 - c. Connect the BNC end of the test cable to the 50 ohm feed-through connector on scope channel 1.
22. **[Manual Prescan] [Scan TR].** Increase TG to 200.
23. Read the peak-to-peak size of the waveform in divisions and then confirm that it is ≥ 35.7 **minor divisions.**
24. Decrease TG to 0 (zero). **[Done].**
25. **If the system has a (U)CERD:** Remove the cannon-to-BNC test cable from the top-most plug of J109 on the (U)CERD and reconnect the cable that was previously removed from J109.
27. **If the system has an RRF:** Remove the SMB-female to BNC-male test cable from the Master Exciter Output J8 connector on the Multiplexer Board and reconnect the cable that was previously removed.
28. **If the system has an RRF:** Trace the cable connecting to the Master Exciter Input J10 on the Multiplexer Board over to the J6 RF Out on the Master Exciter Board.
 - a. Remove the cable at J6 RF Out on the Master Exciter Board and connect the SMB-female to BNC-male test cable to this port.
 - b. Repeat steps 22 – 24 above.
 - c. Remove the SMB-female to BNC-male test cable from the J6 RF Out connector and reconnect the cable that was previously removed.
29. Only consider replacing the Exciter or RRF Master Exciter Boards if the output is below the specification **AND** the body RF output power cannot be adjusted to meet the 16 kW specification. Otherwise, check the RF cabling, connectors or any hardware in between.

APPENDIX F — DUMMY LOAD AND CABLES CALIBRATION

Description - This procedure provides directions for determining the true loss attributable to the dummy load and cables used when measuring the RF output power with a wattmeter or oscilloscope of 100 MHz bandwidth or greater. It is necessary to know and account for the actual loss these components contribute in order to accurately measure RF power. This procedure is not needed if using the RF Power Measurement Kit. The RF Power Measurement Kit has already been calibrated so that this loss is known and accounted for.



WHEN MEASURING HIGH FREQUENCY (IN THIS CASE, 42 OR 64 MHZ) ON ANY 100 MHZ SCOPE, THERE MAY BE AN ERROR DUE TO THE BANDWIDTH LIMITATIONS OF THE SCOPE. (IF NECESSARY, READ 'SCOPE TRIVIA' AT THE END OF THIS PROCEDURE.)

F-1 Overview

Test cables long enough to reach the cables, connectors and dummy load to be tested are connected between the Exciter RF Output (System Cabinet J1) and Receiver Body Input (System Cabinet J2). Receiver gains (R1 & R2) and transmit gain (TG) are set for near full scale reading on the power spectrum during prescan calibration. A reference scan is taken and stored in a raw file. The Attenuation Test Tool is used to calculate the baseline factor from the reference scan for the test cables (i.e., there is some loss from the test cables).

The dummy load and/or cable(s) to be tested are next inserted in series with the test cables and another scan is taken. Again, the Attenuation Test Tool is used to determine the "Magnitude Squared Attenuation Factor" (i.e., how much has the test signal been attenuated?). This attenuation factor is used in the RF power calibration process to accurately calculate the RF power level.

Note

If any problems are encountered during the following procedure, always start over at the beginning and re-do the reference scan. Then you may add, as directed in this procedure, any type of attenuation hardware you might have reason to test.

F-2 Tools and Instruments Required

See Table F2-1

TABLE F2-1
REQUIRED TOOLS AND INSTRUMENTS

Item	Description	Part Number	Qty.
1	50-ohm dummy load, 200 watt, 30 dB attenuator - Bird Model 8322 (or equivalent).	46-317724P14	1
2	RF Test Cables Kit	46-255816G1	1

F-3 Initial Setup

1. Disable the (U)TNS.
 - a. **If system is equipped with a TNS** then bypass the body TNS module with a female-BNC to female-BNC adapter.




Failure to bypass the body TNS out of the circuit may result no signal being received by the system. Merely disabling the TNS by moving the Enable/Disable toggle switch down to the Disable position instead of bypassing it out of the circuit may work, however, the TNS processor can override the switch.

- i. Locate the body TNS on the TNS module. This is one of the two shiny metal TNS boxes affixed to the main TNS module assembly that is mounted farthest from the rear of the cabinet. Multicoil systems have an additional piggyback board affixed with 3 extra TNS boxes that mounts over the top of the main TNS module assembly.
 - ii. Bypass the body TNS (A24 A1 A1) from the circuit by disconnecting the small coax cables A24 A1 A1 J12 (signal output) and A24 A1 A1 J2 (signal input) from the TNS and connecting both together using a female-BNC to female-BNC adapter (also known as a bullet adapter 46-220427P3).
- b. **If the system is equipped with a UTNS (MGD hardware only):** It is not currently necessary to disable the UTNS. No action is required.

F-3 Initial Setup (Continued)

2. Reconfigure test hardware as shown in Illustration F3-1.

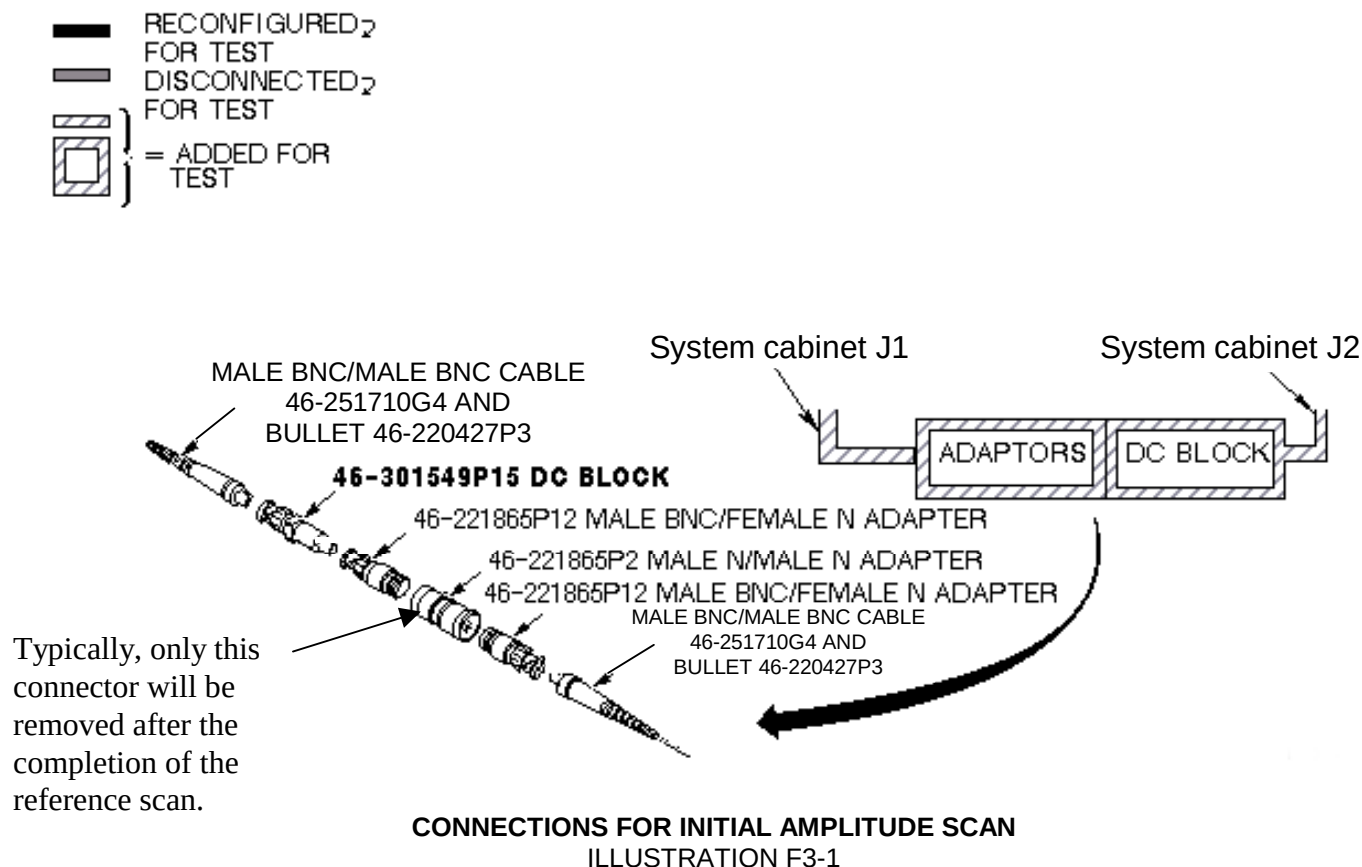

Note

Adapters not shown in Illustration F3-1 can be added, if necessary, from the RF Cables Kit. Usage of the inline DC Block as shown in Illustration F3-1 is mandatory.

3. Disconnect the existing cables at the Systems cabinet I/F panel J1 (RF Out) and J2 (Receiver Body Input) and set them aside.
4. Connect the assembled test cables and adapters between J1 (RF Out) and J2 (Receiver Body Input) on the System Cabinet interface.

F-3 Initial Setup (Continued)

5. At the operator work space, prepare the system for a Dummy Load scan using the procedure, see below.
 - a. Click on **[New Pt]**
Id: **geservice**
Name: **dummy load**
Weight (Lb): **111**
Set Patient Protocols to **Service**.
 - b. In the Protocol field, type **o.18.1** (o=Other, 1=series number) to load the protocol.
 - c. Set a landmark if necessary, then **Save Series**.
 - d. With the right mouse key, select **[Research Operations]**, then select **[Display CVs]**.
Set value of CV **calmode** to **2** (trapezoid pulse).
(Caution here. Make sure the previous CV has been cleared before entering the next one. Look at the screen!)
Set value of CV **p2_ramp** to **1** (1 μ sec ramp time).
Set value of CV **t2** to **50000** (50 msec tr).
Set value of CV **pismode** to **1** (exc service).
Set value of CV **pmode** to **1** (data collection).
Set value of CV **daqm** to **1** (data in window).
 - e. Select **[Accept]** and then select **[Research Operations]**, then select **[Download]** then select **[Manual Prescan]**.

F-4 Data Collection

1. When in **[Manual Prescan]**, set **R1** to **7**, and **R2** to **14**.
2. Adjust transmit gain (**TG**) to achieve an R1 or R2 (on IP display) of approximately 98%, without going over.
3. Select **[Done]**.
4. Select **[Scan]** (Ignore the message: MR signal too large, reduce receiver gain.) (Note: on the LX systems tested, the scan time starts at 13 seconds, counts down to 7 seconds, then ends. This is normal and is not cause for alarm.)
5. From the MR Tools desktop, select **[Cals/checks]** and then **[Attenuation Test]**.
6. Use **[Atten Test]** tool selection to analyze data, as shown in Table F4-1.

TABLE F4-1
DATA COLLECTION

Output/Prompts	Input/Comments
Last run number used was: XXXX Please enter runfile number (XXXX): Please select Locked / Unlocked file (L,U) (U):.....	<Enter> <Enter> (working)
<pre> ***** ***** Average Max. magnitude Across All Views = aaaaa Average Max. magnitude Squared = bbbbb Average RMS Across All Views = ccccc ***** Do you want to make this run the reference(Y,N)(N):..... </pre>	Y<Enter>
STOP! Do not answer the next question at this point. Continue with step 7 below.	

7. The next step will involve removing only the center male N to male N adapter from the test cables, setting it aside, and adding in the dummy load and cables that need to be characterized.



Note

This assumes that the item to be characterized has N female connectors at it's input and output. If it has an N connector at the input and a BNC at the output or BNC connectors at both the input and output then an additional adapter(s) will be needed in order to connect it to the test cables. Adding in one or two uncharacterized adapters should not appreciably change the baseline attenuation factor. In this case, if the wattmeter is not being used, the 2.5 ft. cable (RG214, 46-255595G4) can be eliminated from the circuit. It is not needed.

F-4 Data Collection (Continued)

8. Connect your test cables to the opposite ends of either:

- Illustration F4-1 - Dummy Load and Cables
- Illustration F4-2 - Amplifier to Wattmeter Cable.

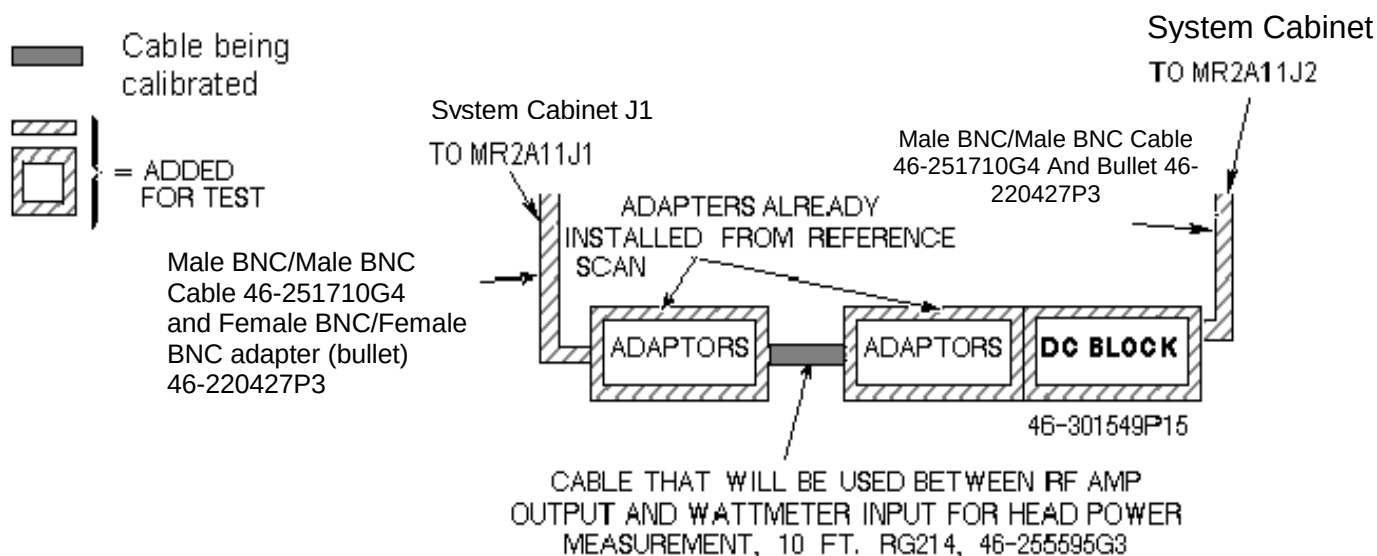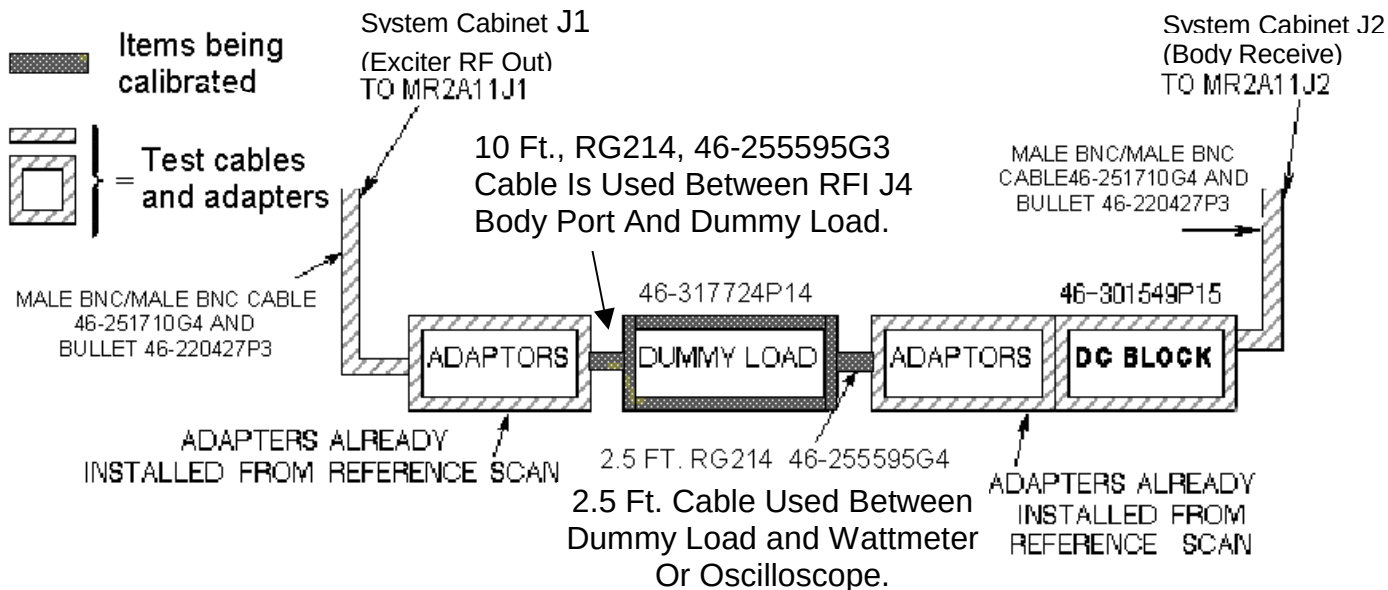

Note

"Bullet" RF connector referred to in illustrations F4-1 and F4-2 is a female BNC/female BNC adapter, 46-220427P3.

F-4 Data Collection (Continued)

1. Select the scanning icon again to activate the scanning screen.
10. Select **[Scan]**.
11. When the scan is completed, re-select the tools icon again and begin the Analysis in section F-5. It will be necessary to toggle between the Scan and Toolbelt icons if multiple passes are done.

F-5 Analysis

See Table F5-1.

TABLE F5-1
ANALYSIS

Output/Prompts	Input/Comments
Do you want to compute Gain or Attenuation Ratio(G,A)[G]:	A<Enter> IMPORTANT: Answer after the scan is done.
Please do the next scan...	
Press <enter> key to continue, s or q to quit []:	<Enter>
Last run number used was: XXXX	
Please enter runfile number	<Enter>
[XXXX]:.....	
Please select Locked / Unlocked file (L,U)	<Enter> (Working)
[U]:.....	

Average Max. Magnitude Across All Views = aaaaa	
Average Max. Magnitude Squared = bbbbb	
Average RMS Across All Views = ccccc	
Magnitude Attenuation Factor = xxxxx	
Magnitude Squared Attenuation Factor = yyyyy	<=====Record in "Value" column in Table F5-2.
RMS Attenuation Factor = zzzzz	
Do you want to make this run the reference(Y,N)(N):.....	N<Enter>

1. Record the "Magnitude Squared Attenuation Factor" number and record it in the appropriate Value box in Table F5-2.

TABLE F5-2
ATTENUATION FACTORS

MODE	CALIBRATED HARDWARE	PART NUMBER(S)	VALUE	NOMINAL VALUES
BODY OR HEAD	DUMMY LOAD + CABLES ATTEN FACTOR	46-255595G3 46-317724P14 46-255595G4		1000 TO 1200
HEAD	AMP TO WATTMETER ATTEN FACTOR	46-255595G3		1.0 TO 1.2

F-6 Calculation of RF power

Peak voltage should be used in the calculation in order to get an accurate result. It can be converted to power using the following formula as long as certain factors are known and accounted for. The scope correction factor **MUST** be known. So must the *actual* total loss attributed to anything that connects the measuring device to the source. This often includes the accumulated loss associated with the dummy load and any interconnecting cables. Table F6-1 shows the calculation of power if all the attenuating devices in the measurement circuit exhibited perfect loss; that is, the devices added no more or less loss than what they were designed to provide. Table F6-2 shows the same calculation of power but accounts for the measurement-circuit loss values in deriving the true power. Note that the loss has a significant impact on the calculated power.

WARNING!

THE SCOPE CORRECTION FACTOR MUST BE KNOWN FOR THE FORMULAE SHOWN IN TABLES F6-1 AND F6-2. IF IT IS NOT KNOWN, DO NOT USE THIS METHOD. GROSSLY INACCURATE MEASUREMENTS AND POSSIBLE SYSTEM DAMAGE WILL RESULT.

$$\frac{\left(\frac{V_{\text{peak}}}{\text{scope correction factor}} \right)^2}{2 \times Z} \times \text{dummy load and cables attenuation} = \text{RF Power}$$

where "X" and "X" in the above formula signifies multiplication

Assume $Z = 50 \, \Omega$, $V_{\text{peak}} = 40.0$, scope correction factor = 1.00 (no loss), dummy load and cables atten. = 1000, (dummy load and cables are all ideal)

$$\frac{\left(\frac{40.0V_p}{1.00} \right)^2}{100} \times 1000 = 16000 \text{ Watts} = 16 \text{ kW}$$

This result assumes a ***theoretically perfect*** situation in which there is no loss. These situations, in common practice, rarely exist!

TABLE F6-1
RF POWER CALCULATION WITH NO LOSS

F-6 Calculation of RF power (Continued)

Now, consider the "real life" type situation in Table F6-2 in which the loss is considered:

$$\frac{\left(\frac{V_{\text{peak}}}{\text{scope correction factor}} \right)^2}{2 \times Z} \times \text{dummy load and cables attenuation} = \text{RF Power}$$

where "X" and "X" in the above formula signifies multiplication

Assume $Z = 50 \, \Omega$, $V_{\text{peak}} = 34.61$, scope correction factor = 0.88, dummy load and cables atten. = 1028

$$\frac{\left(\frac{34.61V_p}{0.88} \right)^2}{100} \times 1028 = 15901 \text{ Watts (very close to 16kW.)}$$

Accounting for the loss resulted in an accurate answer. 15899 Watts is as close as we can hope to get to 16000 Watts without using the RF Power Measurement Kit. Note that if none of the loss had been accounted for the error could have been **3587 Watts or 22.6%**. As a result, the observer would attempt to adjust the RF power far above the 16000 Watt limit!

TABLE F6-2
RF POWER CALCULATION WITH LOSS CONSIDERED

F-7 System Restoration

1. Reconnect original cables to the System Cabinet I/F J1 and J2.
2. Enable the (U)TNS.
 - a. **If system is equipped with a TNS** then remove the female-BNC to female-BNC adapter and reconnect RF cable A24 A1 A1 J12 (signal output) to the bottom of the body TNS module and RF cable A24 A1 A1 J2 (signal input) to the top of the body TNS module.
 - b. **If the system is equipped with a UTNS:** It is not currently necessary to disable the UTNS. No action is required.
3. Perform one satisfactory head or body scan.

F-8 Wattmeter Trivia

Why are measurements between a oscilloscope and a Bird Wattmeter found to be different? Repeated GEMS Education Center lab groups that performed the Dummy Load and Cables attenuation process correctly, and obtained an accurate attenuation value, were able to accurately measure RF power as well as with a 400 MHZ scope! Remember, for wattmeter measurements however, the meter reading must be multiplied by the atten factor. Leaving this little factor out is the reason why the above 'note' was made. In body mode the wattmeter is placed after the dummy load. Its atten value must be used to calculate the correct RF value. If it is NOT used, then it should be expected there will be an error of 10% or more!

In head mode, the wattmeter is placed before the dummy load. In this case only the 'cable' factor is used to correct the wattmeter reading. The 'cable' factor is the other atten test you did for the 10 foot cable only. Typically that value is between 1.0 to 1.2. Doesn't sound like much???? Take your calculator out for a moment. Take 2,000 watts (displayed reading) and multiply by an atten factor of 1.0 ; the results are obvious, 2,000 watts. Now take the same displayed 2,000 and multiply it by an atten factor of 1.2 = 2,400 watts! Hopefully this example will show the importance of the attenuation values. You just calculated a 20% error if you were wrong because you didn't do atten test on the 10 foot cable!

F-9 Scope Trivia

A basic 'rule of thumb' in the scope industry is that 'to accurately measure a high frequency signal, the bandwidth of the scope should be at least 5X the frequency being measured'. That means that for a 64 MHZ signal the scope bandwidth would have to be at least 320 MHZ. Along the same line, a 42 MHZ signal would require a bandwidth of at least 210 MHZ. This is why the 400 MHZ scope was chosen several years ago for MR use. Knowing this should alert you as to what might have to be done in order to get correct scope readings from whatever scope you carry today. Who knows what bandwidth scope you'll be working with tomorrow? Will it be of sufficient bandwidth to give accurate measurements? What can you do if it is not covered by the rule of thumb? Read on!

Experimentation with multiple scopes have shown that an error is typically somewhere between 0 and 16%. That's right, some 100 MHZ scopes have no loss! Others may have as much as 15 or 16%. Additionally, it has been found that it is not necessarily the whole scope that has this loss, but the individual channels typically exhibit differences. In other words, channel one might have a 3% loss when channel two could lose 15%. There is no fixed number for all 100 MHZ scopes!

How do you get around these losses? Well, there are two ways. First you could compare a RF waveform on your scope, side by side, with a 400 MHZ scope. Use the same cables etc. and just move the cables from scope to scope. You can then calculate the loss for each channel of the 100 MHZ scope vs. the 400 MHZ scope. Keep the scope intensity low for best accuracy. If and when possible, determine the losses for a 1.5T and 1.0T system. That's for 64 and 42 MHZ. Point of note: there are various losses at 42 MHZ as well as with 64 MHZ. And, typically they differ! Your scope may have a 0 to 15% rolloff at 64 MHZ and 0 to 12% at 42 MHZ. All scopes differ!

Secondly, the next time you have your scope calibrated, just attach a tag asking for the % loss of signal at 42 and 64 MHZ for each channel. You should also state that you want this to be done at the 5 and 2 volts/div setting. This will cover both body and head measurements. At two different re-calibration vendors we use, there are no extra charges incurred for this request.

As mentioned above, the 100 MHZ scope has been shown to be as accurate as a 400 MHZ scopes if the scope rolloff factor is known. Additionally, each channel has it's own unique loss factor.

The biggest error found in GEMS Education Center labs when measuring RF power has proven to be the FE using "typical" dummy load and cables attenuation factors. Basically there are no "typical" factors! All loads and cables are somewhat different. Several may appear to be close in value, but in most cases they are different. If you use a "typical" atten factor in your calculations, then "typically" you will be in error.


Through experimentation at the GEMS Education Center, it has been determined that if you know the attenuation factor for a particular dummy load and cables, those values hold over several years. They do NOT tend to change. This means that if you could do atten test on a dummy load and cables, and keep them as a kit, you don't have to recalibrate them every time you use them. Wow! And possibly a little cheaper than the Power Measurement kit.


F-9 Scope Trivia (Continued)


If you must, or intend to use a scope to measure body RF power, use atten test on the dummy load and cables.

Knowing the correct attenuation values and the correct scope rolloff values will provide you with measurements as accurate as a 400 MHZ scope and / or the GE Power Measurements kit.

APPENDIX G — GENERAL TROUBLESHOOTING

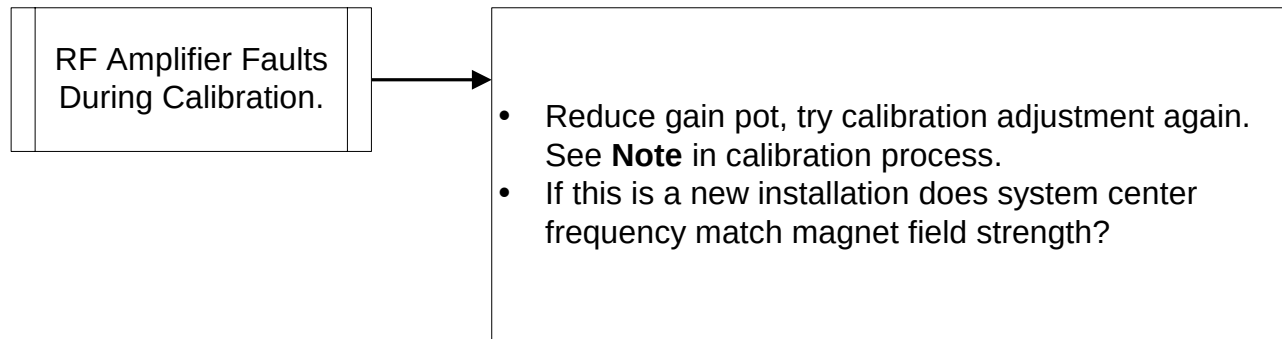

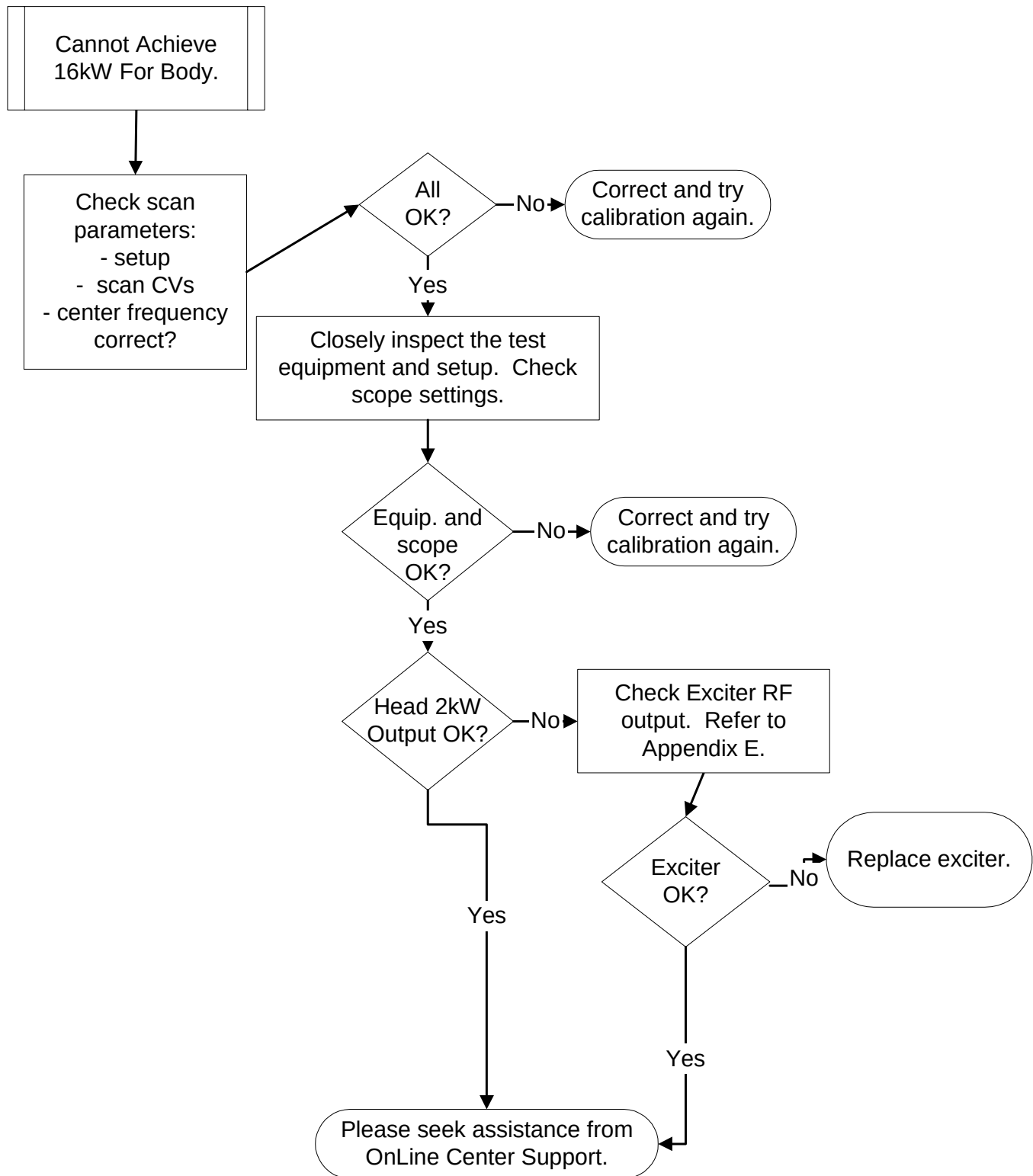

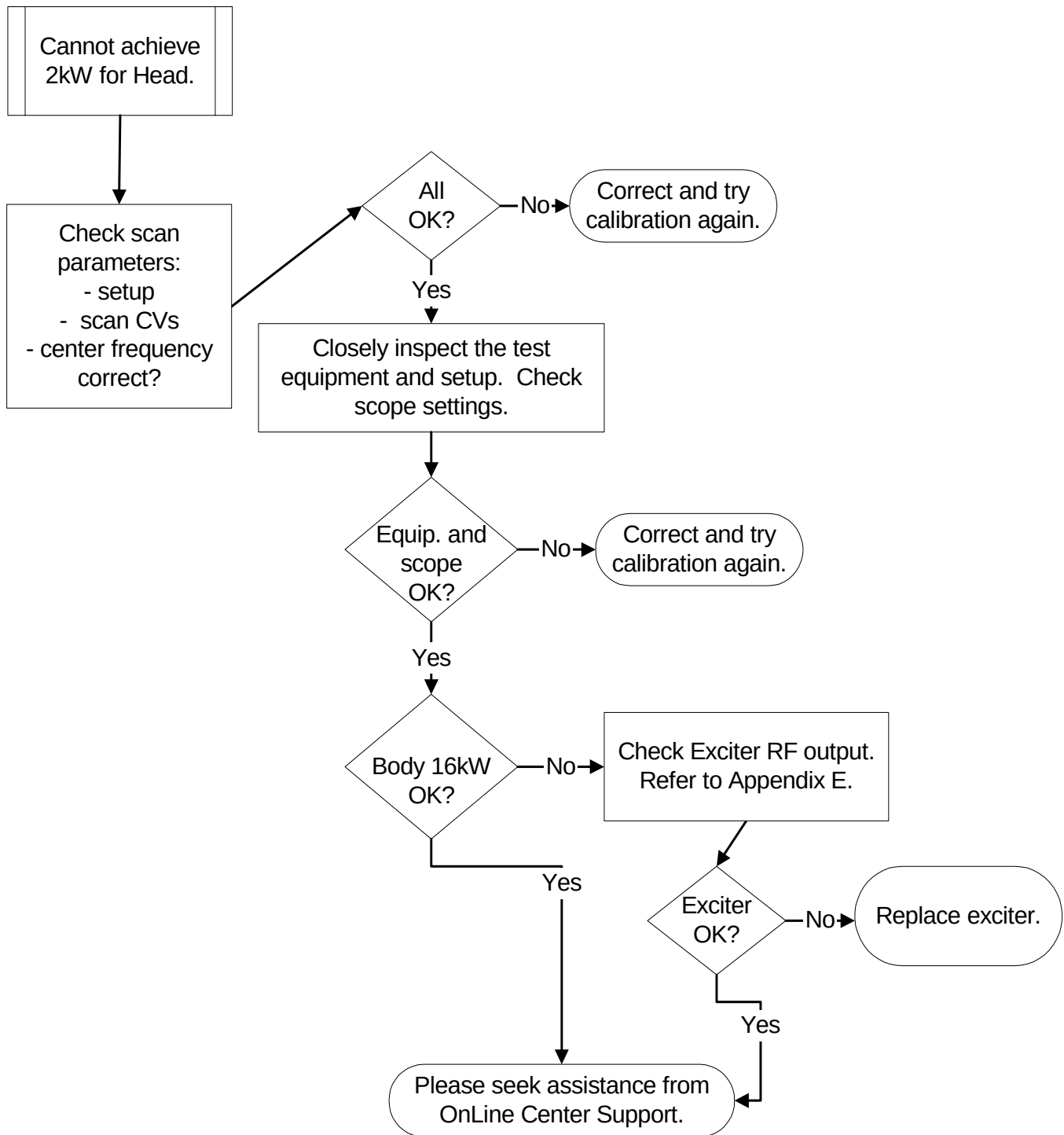

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	October 23, 2002	K. Keshena	New SRFD2 calibration document.
1	November 11, 2002	D. Thome	Added Item E "Prepare to Scan" to Tables 2-4 and 3-2. Modified document to accommodate MGD Hardware.