

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1 - OVERVIEW.....	2
2 - TOOLS AND INSTRUMENTS REQUIRED.....	4
3 - INITIAL CONDITIONS.....	5
3-1 Pre-requisite Status.....	5
4 - INITIAL HARDWARE PREPARATION	6
5 - HEAD HARDWARE SETUP AND CHECKS.....	11
5-1 Head Hardware Setup.....	11
5-2 Converting Units If Using The RF Power Measurement Kit ONLY	12
5-3 Head Scan Setup	16
5-4 Head Power Monitor Functional Check Process.....	16
6 - BODY HARDWARE SETUP AND CHECKS	22
6-1 Body Hardware Setup	22
6-2 Converting Units If Using The RF Power Measurement Kit ONLY	23
6-3 Body Scan Setup.....	27
6-4 Body Power Monitor Functional Check Process	27
7 - AMPLIFIER SHUTDOWN VERIFICATION	33
8 - SYSTEM RESTORATION	36
9 - DECLARATION FORM PREPARATION	37
APPENDIX A — SIGNA POWER MONITOR FUNCTIONAL TEST DECLARATION FORM	38
APPENDIX B — ALTERNATE EQUIPMENT SETUPS & NON-SERVICE PROTOCOLS.....	39
APPENDIX C — SOFTWARE AND FORMULAE FOR MISCELLANIOUS RF CONVERSIONS	42
C-1 RF Conversion Software.....	42
C-1-1 RF Power Calculator.....	42
C-1-2 Power Calculator Tool.....	42
C-2 RF Conversion Formulae.....	42
APPENDIX D — DUMMY LOAD AND CABLES CALIBRATION.....	44
D-1 Overview.....	44
D-2 Tools and Instruments Required	45
D-3 Initial Setup	45
D-4 Data Collection	48
D-5 Analysis.....	51
D-6 Calculation of RF power	52
D-7 System Restoration	53
D-8 Wattmeter Trivia	54
D-9 Scope Trivia.....	55
REVISION HISTORY	57

1 - OVERVIEW

This section applies to systems installed with the 1.0T RF/PDU Cabinet (contains the RFI, SSM, 2 Teal RF amplifiers, and PDU).

Use this procedure to ensure that the power monitors are fully functional, including redundant trip capability. The tests described use a special PSD (pulse sequence data base) to provide excitation for the RF amplifier. The PSD allows the service engineer to selectively force peak power, pulse width, and duty cycle beyond allowable limits. Select set points for testing by modifying Control Variables that override the values normally calculated and downloaded by the PSD.

Power Monitor firmware operation is changed to bypass/test mode by moving “monitor” switches (A and B) on the front of the SSM. This action permits set-point verification without tripping the RF amplifier, thereby eliminating the time consuming necessity of resetting the amplifier after each trip. It is necessary to measure the RF output power and know at what level each monitor tripped.

Three methods are available and listed in the order of preference for measuring head and body RF output power:

- **RF Power Measurement Kit** - Easy to use, preferred method of power measurement.
- **Bird wattmeter** - Must account for cable and load loss to determine actual power.
- **Oscilloscope** - Must account for cable and dummy load loss to accurately calculate the power from the observed peak voltage. This method works but the process used to derive the power is somewhat complex and lacks the accuracy of the previous two methods. It is, for this reason, not recommended.

It is *strongly recommended* that the RF Power Measurement Kit be used to obtain an accurate measurement of the RF output power. If it is not possible to obtain an RF Power Measurement Kit then one of the two other methods can be used provided that:

- The Dummy Load and Cables Calibration procedure in Appendix D has been performed and the individual loss values are known.


WARNING!

FAILURE TO ACCOUNT FOR THE ACTUAL LOSSES IN THE DUMMY LOAD AND CABLES WHEN NOT USING THE RF POWER MEASUREMENT KIT TO MEASURE RF POWER CAN RESULT IN SIGNIFICANT MEASUREMENT ERRORS.

After checking the set points, the actual shutdown of the amplifier is verified: once for monitor A, once for monitor B, and once for the HV relay shutdown. The actual shutdown of the amplifier is verified twice for monitor A (to verify both relay and logic shutdowns), and once for monitor B (connected to amplifier via monitor A). This also confirms that the bypass/test “monitor” switches (A and B) are restored prior to returning the system to normal operation.

A Calibration Flowchart has been provided in Illustration 1-1.

1 – OVERVIEW (CONTINUED)

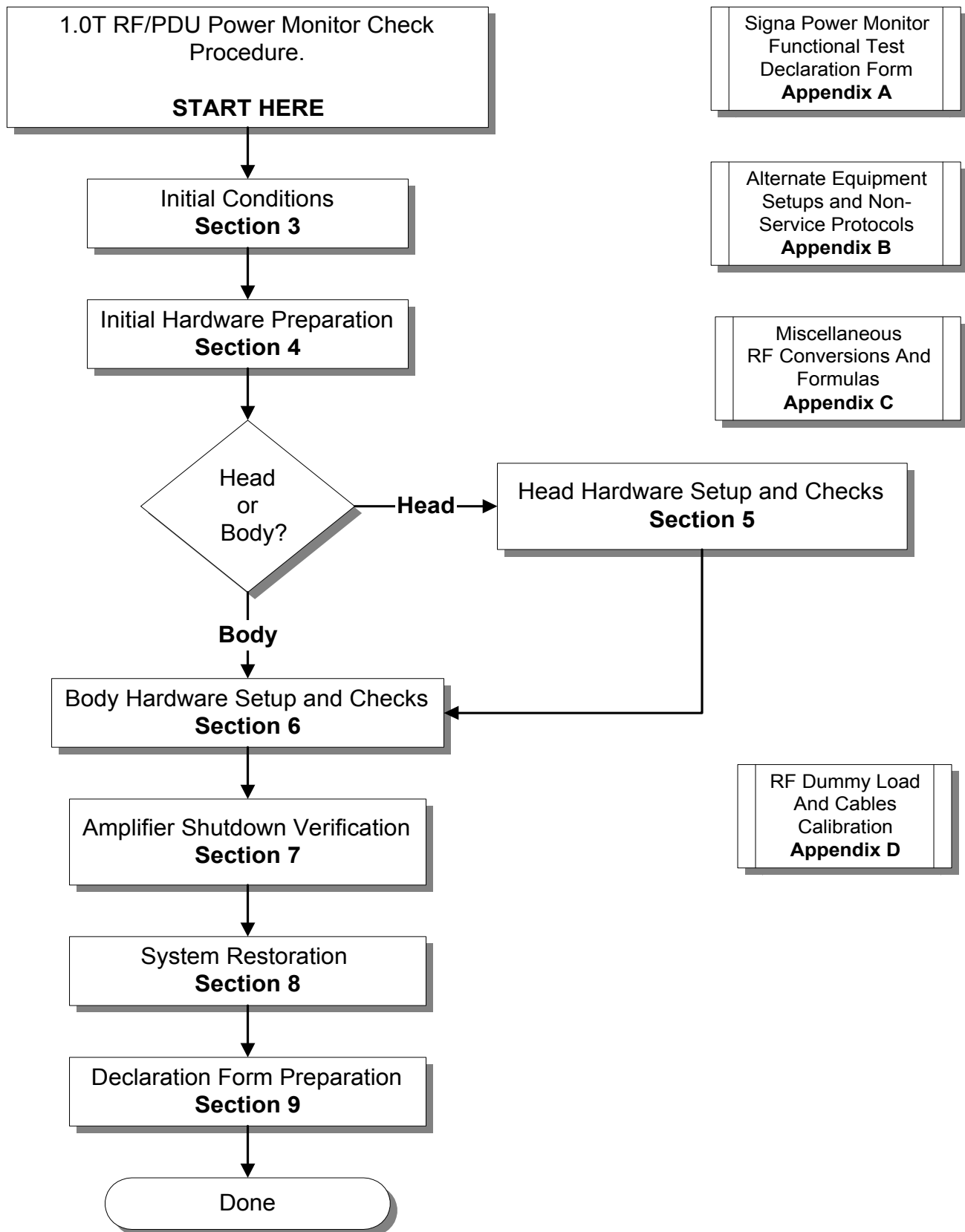

POWER MONITOR CHECK FLOWCHART
ILLUSTRATION 1-1

2 - TOOLS AND INSTRUMENTS REQUIRED

A Laptop and Serial Port I/F cable (vendor part # 510074) will be required if performing the proprietary procedure. Table 2-1 describes what tools are necessary if the RF Power Measurement Kit is NOT used. Table 2-2 describes all the necessary tools needed if the RF Power Measurement Kit is used.

The use of the RF Power Measurement Kit to complete the RF Power Monitor Check is highly recommended.

TABLE 2-1
TOOLS AND INSTRUMENTS REQUIRED WITHOUT RF POWER MEASUREMENT KIT

Item	Description	Part Number
1.	RF Test Cables Kit	46-255816G1
2.	50 ohm, 200 Watt, 30 dB attenuator - Bird Model 8322	46-255837P10
3.	Oscilloscope	46-183029P61
4.	Laptop serial cable (optional)	2124497-46
4.	Wattmeter and appropriate elements (optional)	Not supplied

TABLE 2-2
TOOLS AND INSTRUMENTS REQUIRED WITH RF POWER MEASUREMENT KIT

Item	Description	Part Number
1.	RF Power Measurement Kit	46-317724G1 or G2
2.	50 ohm, 200 Watt, 30 dB attenuator - Bird Model 8322 NOTE: Only required with G1 Kit	46-255837P10
3.	Oscilloscope	46-183029P61
4.	Laptop serial cable (optional)	2124497-46

3 - INITIAL CONDITIONS

Perform the prerequisite calibration procedures in the sequence listed in Table 3-1.

TABLE 3-1
PRE — REQUISITE PROCEDURES

Step	Procedure Title	When Required
1	Dummy Load & Cables Calibration (Do ONLY if not using the RF Power Measurement Kit)	Calibration required whenever using the dummy load and cables. This value will be frequency specific and can change over time.
2	RF/PDU Max Power RF Out Setup and Calibration	At initial installation or if the RFI, UCERD Exciter, or RF Amps have been repaired or replaced.



PREVENT EQUIPMENT DAMAGE. ENSURE THAT THE RF AMPLIFIER HAS A DUMMY LOAD CONNECTED TO ITS OUTPUT BEFORE INITIATING TESTS. IF NOT, FAULTS MAY OCCUR AND RF AMPLIFIER COMPONENTS COULD BE PERMANENTLY DAMAGED.

3-1 Pre-requisite Status

Pre-requisite — Signa must be fully functional and properly calibrated.

Pre-requisite — Software must be running with the system capable of normal imaging.


4 - INITIAL HARDWARE PREPARATION

1. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the **DANGER** message on this page.




PROPERTY DAMAGE! TO PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, REMOVE ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

2. Remove the front and rear covers from the RF/PDU Cabinet.
3. See Illustration 4-1. On the front of the SSM place the:
 - 2 (two) power MONITOR switches (A and B) to the middle BYPASS position.
 - 3 (three) DRV FLT switches to the bottom DIS (disable faults) position.

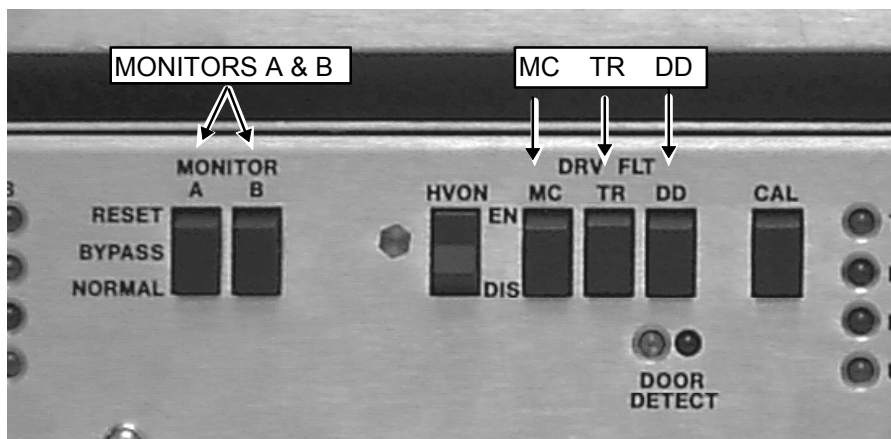

SSM FRONT PANEL SWITCHES DISABLED
ILLUSTRATION 4-1

4. Remove the test window from the Test Switches and Test LEDS on the front of the SSM. See Illustration 4-2

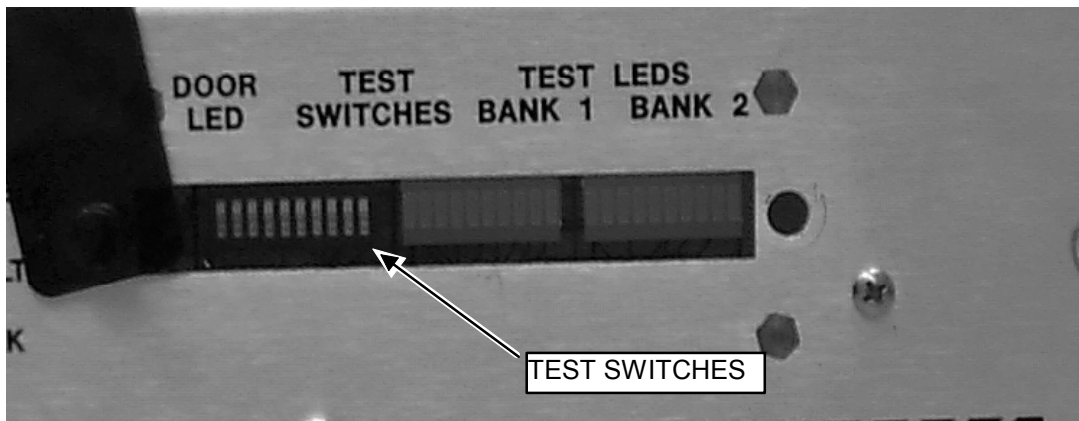

INSIDE THE TEST WINDOW
ILLUSTRATION 4-2

4 - INITIAL HARDWARE PREPARATION (CONTINUED)

Note

Order of DIP switch placement is important otherwise the system will not recognize changes.

5. Place the first four switches in the order and positions shown in Table 4-1 and leave in these positions for the duration of the test.

TABLE 4-1
SWITCH POSITIONS

SWITCH NUMBER	POSITION
2	UP
3	DOWN
4	UP
1	UP

6. The TEST LED will illuminate. This is the yellow LED to the right of the CAL switch.

4 - INITIAL HARDWARE PREPARATION (CONTINUED)


7. Refer to Table 4-2 for the alternate proprietary procedure hardware configuration.

Note

The alternate proprietary procedure utilizes the field laptop and the MONS1.EXE software to detect any faults from RF power monitors A and B. A long (length depends on the distance between the RF/PDU Cabinet and the OW Console at your site) 25 pin female to 25 pin male serial cable (NOT NULL MODEM) can be used with cable 2124497-46 described below in Table 4-2 to remotely monitor and reset power monitor faults from where prescanning changes are being made at the OW Console. While the cable length probably shouldn't exceed 50 feet (15.24 meters) in order for reliable remote monitoring to take place, longer cable lengths (especially if the cable in use is shielded) may work. Alternatively, if the 2124497-6 cable is unavailable, the same locally procured long cable with a null modem adapter affixed to one end and a 25 pin female to 9 pin female adapter connected to the other end could also be used. See Table 4-2 below. Refer to **SECTION 3 - MONS1.EXE Revision 2.2 in LAPTOP CONFIG AND TROUBLESHOOTING SOFTWARE** (RC3TSC5.DOC) on the 8.X Service Methods CD-ROM for more detailed information concerning the use of MONS1.EXE software.

TABLE 4-2
ALTERNATE PROPRIETARY PROCEDURE USING LAPTOP

Connect the laptop serial port to J46 on the rear of the SSM, using cable part number 2124497-46 (Vendor part # 510074, supplied with every cabinet). See Illustration 4-3. If this cable is not available, use a suitable null modem and 9-25 pin (female on both ends) cable adapter.

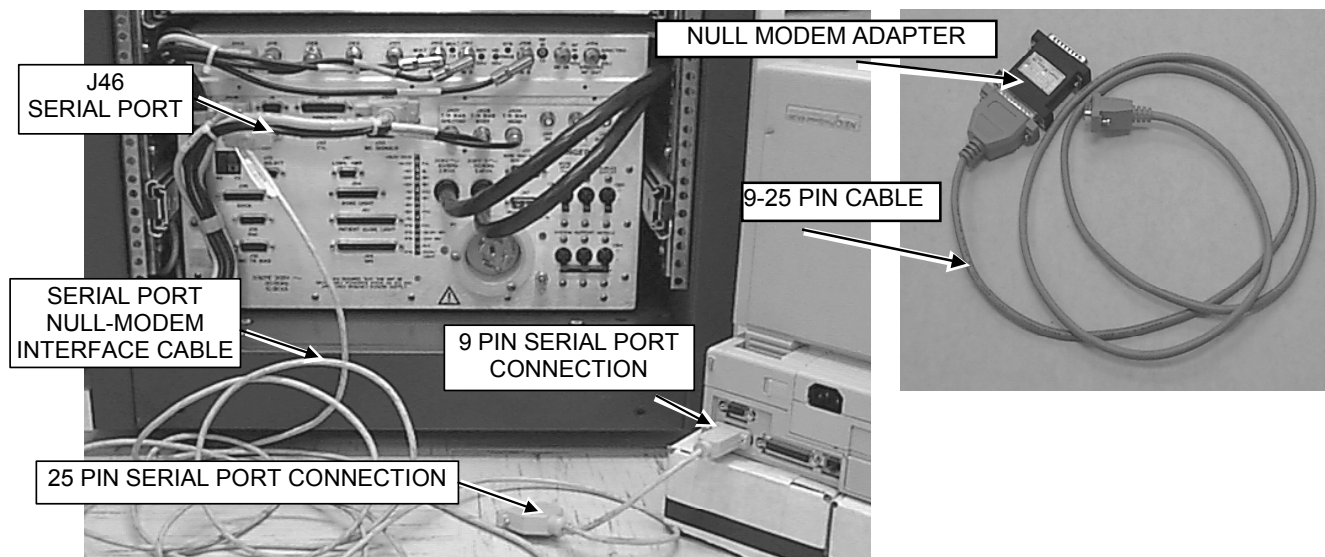

4 - INITIAL HARDWARE PREPARATION (CONTINUED)

8. If using the alternate proprietary procedure, see Table 4-3 for instructions for running the MONS1.EXE software.

Note

Only version 2.2 of MONS1.EXE is compatible with all MR system types. This is the version currently available from the Signa LX Service Methods CDROM, revision 12 and later.

TABLE 4-3
ALTERNATE PROPRIETARY PROCEDURE — LAPTOP SOFTWARE

Start the software diagnostic program MONS1.EXE on the laptop.

Note

The MONS1.EXE program, if not loaded on the service laptop, can be loaded from the Signa LX MR Service Methods CD-ROM 2160623-1. Insert the CD-ROM into the drive. Under Windows 95/98 use the cursor to select "**Start**" and then "**Run**" and then "**Browse**". From the Browse window, select the CDROM drive from the down arrow to the right of the "Look in:" entry box. Select the **mrtools** directory, then highlight the **Setup.exe** file and select "**Open**" at the bottom of the window. Next, select "**OK**" from the bottom of the Run window.

Follow the setup program instructions to install the tool onto the laptop. The MONS1.EXE must be run in MS-DOS mode and not from a Windows shell to DOS. Once the installation is complete, from Windows 95/98 perform the following steps:

- 1) **[Start], [Shut Down...], Restart the computer in MS-DOS mode?, [Yes]**

Note

It may take up to 1.5 minutes for some Dell laptops to restart in MS-DOS mode.

- 2) At the C:\Windows> prompt, type **cd \cclass\rfertools<ENTER>**
3) At the C:\CCCLASS\ERBTEC> prompt type **mons1.exe<ENTER>**

9. Once the **mons1.exe** program has started then select **N** (next screen) and then **T** (test mode) to reach a screen like that shown in Illustration 4-3.

	MONITOR A	MONITOR B
A Fault is	NOT DETECTED	NOT DETECTED
Press 'esc' to return to previous screen, 'C' to Continue after Fault		

'T' TEST MODE SCREEN (NO FAULT)
ILLUSTRATION 4-3

The screen shown in Illustration 4-3 will be used to remotely and individually monitor Power Monitors A and B for faults.

4 - INITIAL HARDWARE PREPARATION (CONTINUED)

When a fault occurs the word **DETECTED** will be displayed below the power monitor that registered the fault. Illustration 4-4 shows a situation where both power monitors registered a fault.

	MONITOR A	MONITOR B
A Fault is	DETECTED	DETECTED
Press 'esc' to return to previous screen, 'C' to Continue after Fault		

'T' TEST MODE SCREEN (WITH FAULT)

ILLUSTRATION 4-4

The **C** key on the laptop permits the FE to reset the status of Monitors A and B displayed on the Test Mode Screen to "Not Detected" after the FE has manually reset the Power Monitor A and B Reset Switches on the front of the SSM.

10. **If you are not using the RF Power Measurement Kit** then complete the Dummy Load and Cables Calibration procedure in Appendix D.



When appropriate use the RF Power Measurement Kit, Wattmeter, or Oscilloscope to measure the RF output power and validate that the power monitors are faulting or not faulting at the proper RF power levels. While the Main Screen of the MONS1.EXE software also provides a power measurement readout display for each power monitor, the integrity of each power monitor should be validated independently using the RF Power Measurement Kit, Wattmeter, or Oscilloscope.

5 - HEAD HARDWARE SETUP AND CHECKS

Tests in this section allow the service engineer to selectively force head mode peak power, pulse width, and duty cycle beyond allowable limits.

5-1 Head Hardware Setup



PROPERTY DAMAGE! PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, BY REMOVING ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

1. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the two **DANGER** messages on this page.



PERSONAL INJURY! PREVENT POSSIBLE RF BURNS WHEN DISCONNECTING HELIAX CABLES FROM J3 OR J4 ON THE RFI BY VERIFYING THAT THE SYSTEM IS NOT MANUALLY PRESCANNING OR SCANNING. VERIFY THAT THE SCAN DESKTOP ICON DISPLAYS THE "IDLE" MESSAGE.



2. **If using the RF Power Measurement Kit** then refer to the RF Power Measurement Kit laminated card set.
 - a. Look in the upper, right corner of each card and find the card labeled **58** (1.0T Head RF Output).



The head RF output connection is no longer to the non-existent EFB unit, as the reference cards in some of the older kits show, but instead to the J3 output on the rear of the RFI.

- b. Configure the system as shown in the illustration on the card.
 - c. Confirm that the rotary attenuator is set to the correct position indicated on the card.
 - d. Set the oscilloscope Volts/Div setting of the channel being used to 50 mV/Div.
3. **If using the wattmeter or scope** (NOT the RF Power Measurement Kit) to measure power then refer to Appendix B (Alternate Equipment Setup) for the proper system head configuration.

5-2 Converting Units If Using The RF Power Measurement Kit ONLY

Perform the conversions described in this section only if using the RF Power Measurement Kit. The measurements read from the scope face will be in Volts peak (Vp). The voltage read from an accurate oscilloscope will vary from kit to kit due to variances in cable, attenuator, and directional coupler losses. It is probably easier for most people to determine the expected Volts peak values based on the unique losses of the kit at this point than to wait and convert them later while the values are being read from the scope. The more experienced engineer with a good understanding of RF power measurement principles may, however, prefer to convert and determine the values as each is read from the scope face. This engineer can skip this section and instead refer to **APPENDIX C — MISCELLANEOUS RF CONVERSIONS AND FORMULAS** for the various mathematical formulas needed to perform the conversions.

1. Refer to the RF Power Measurement Kit reference card **58** calibration sticker (1.0T Head RF Output).
2. Record the absolute value of the **Attenuation** number (this means record it without the negative sign), listed on the calibration sticker affixed to the front of the reference card, in the specified location on Table 5-1. This is the total attenuation between the scope channel input and the J3 Head RFI RF port that is provided by the RF Power Measurement Kit. Illustration 5-1 below provides an example of what the calibration sticker looks like.



The **Attenuation** data in the example will be different from what is printed on the calibration sticker affixed to the card in your kit. Do not use the **Attenuation** data listed in the example.

1.0T HEAD CARD # 58 CALIBRATION STICKER DATA
A. Card Number 46-317724P32
Nominal: 58.75 dBm
Attenuation: -55.43 dB
Scope Power: 3.32 dBm
Scope Rdg: 0.927 Vpp

SAMPLE CALIBRATION STICKER
ILLUSTRATION 5-1

5-2 Converting Units If Using The RF Power Measurement Kit ONLY (Continued)

TABLE 5-1
POWER MEASUREMENT KIT HEAD ATTENUATION (IN dB)

ATTENUATION	_____dB
--------------------	---------

3. Start the LX system calculator (**NOT the RF Calculator Tool**) by moving the mouse cursor into the background, right mouse-click, and then select the calculator from the top of the Root menu list.
4. Use the LX system calculator to subtract the Attenuation value in Table 5-1 from the number in each of the 8 locations in Table 5-2.
 - a. Record the result in each of the 8 locations in Table 5-2.
 - b. **Be sure to note the correct polarity sign with each answer. This is critical.** See the example column in Table 5-2.


TABLE 5-2
HEAD FUNCTIONAL CHECK LIMITS (IN dBm)

	Example	Min.	Nom.	Max.
Head Peak Power High Limits in dBm	52 – 55.43 = - 3.43	52 – Attenuation = _____	53 – Attenuation = _____	53.8 – Attenuation = _____
Head Peak Power Low Limits in dBm	50 – 55.43 = - 5.43	50 – Attenuation = _____	50.7 – Attenuation = _____	51.34 – Attenuation = _____
Head Pulse Width High Limits in dBm	46.9 – 55.43 = - 8.53	46.9 – Attenuation = _____	N/A	48.45 – Attenuation = _____

5. Convert each of the 3 Head Peak Power High Limit values (**Min.** [Minimum], **Nom.** [Nominal], and **Max.** [Maximum]) expressed in dBm in Table 5-2 to Volts peak (Vp) using one of the three methods listed below and then record the results in the locations provided in the **Verify** column located in the first row in Table 5-4.
 - a. **METHOD 1:** Use the RF Power Calculation Tool provided on the 8.X Service Methods CD-ROM to convert from dBm to Vp.
 - i. Confirm that the 8.X Service Methods CD-ROM is in the laptop CD-ROM drive and then select the [Power Calculator](#) hyper-link to automatically start the tool.


5-2 Covertng Units If Using The RF Power Measurement Kit ONLY (Continued)

- ii. Alternatively, the tool can be started manually from Windows 95/98 by first moving the cursor to the lower left part of the screen and then selecting "**Start**" .

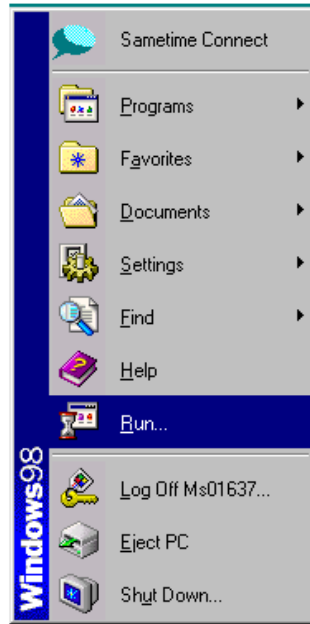

- iii. Next, select "**Run**"

- iv. Finally, type in the file pathname of **E:\rf\power\pwrcalc.htm** as shown:

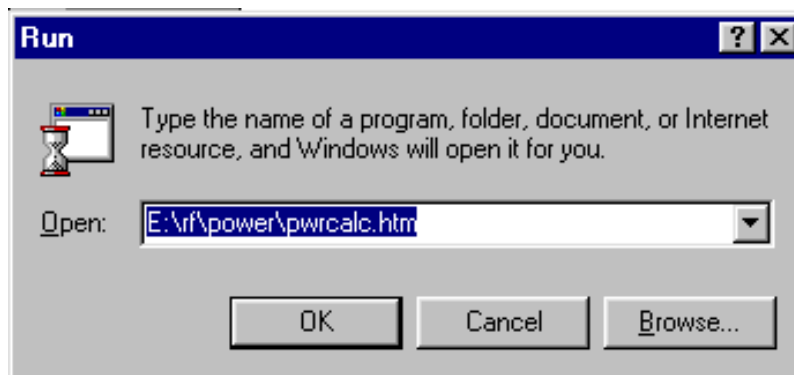

- v. Follow the directions provided in the Tool to perform the dBm to Vp conversion.
- b. **METHOD 2:** Use the LX system calculator to convert from dBm to Vp.
- i. Access the LX system calculator by moving the mouse cursor into the background, right mouse click, and then select the calculator from the Root menu list.

5-2 Coverting Units If Using The RF Power Measurement Kit ONLY (Continued)

- ii. Enter the data into the LX system calculator exactly in the order as shown below:

dBm [=] 20 [=] [INV LOG] [*] 0.316 [=] Vp (←This is the answer)

Example: Using - 3.43 dBm, the result of the first example calculation in Table 5-2, and the LX system calculator to determine the equivalent Vp:

- 3.43 [=] 20 [=] -0.1715 [INV LOG] 0.6738 [*] 0.316 [=] **0.213 Vp**

- c. **METHOD 3:** Use a hand calculator or the scientific calculator available on the field laptop and either of the 2 formulae listed below to convert from dBm to Vp:

Note

“antilog” may also be noted on the laptop calculator as “inv log”.

$$Vp = 10^{\left(\frac{dBm}{20}\right)} \times 0.316 \quad \text{or} \quad Vp = \text{antilog} \left(\frac{dBm}{20} \right) \times 0.316$$

where “X” signifies multiplication

6. Repeat this conversion process for the 3 Head Peak Power Low Limit values in Table 5-2 and record the results in the locations provided in the **Verify** column located in the second row in Table 5-4.
7. Repeat this conversion process for the 2 Head Pulse Width High Limit values in Table 5-2 and record the results in the locations provided in the **Adjustments** column located in the third row in Table 5-4.

5-3 Head Scan Setup

1. At operator workspace, select the scan desktop ICON in the desktop control panel, if you have not already done so.
2. At the operator workspace, prepare the system for a Head Power Monitor scan using the "Service Protocols" procedure in Table 5-3 below. Refer to Appendix B for the non-proprietary protocol.

TABLE 5-3
SCAN PRESCRIPTION - HEAD PM CHECKS

Note: This is the alternate proprietary procedure available for GE use, and to sites with a valid Advanced Service Package Limited License. Refer to Appendix B for the non-proprietary protocol.

- A. **[New Pt]**
Id: **geservice** <ENTER>
Name: **pmc**
Weight (Lb.): **300**
Set Patient Protocols to **Service**.
- B. At front enclosure:
Landmark in the Head area—remove all coils.
press **LANDMARK**.
press **MOVE TO SCAN**.
- C. At Patient Protocols – select **other**.
- D. In the protocol field, type **o.23.1**<ENTER> (o=Other, 23.1 =series) to load the head protocol
OR select **[o.23] [Series 1] [Accept]**.
[OK] (if required).
- E. **[Save Series]**.
- F. **[Prepare to Scan]**.
- G. **[Research Operations]**.
[Setup Params]. Set **TG** to **50**. **[Done]**.

5-4 Head Power Monitor Functional Check Process

1. Refer to Table 5-4. Note that the table consists of 5 columns (**Test**, **Scan Conditions**, **Adjustments**, **Verify**, **Conclusion**) and 10 individual rows of various head power tests.

Note

It may be advantageous at this point to print a hard copy of Table 5-4. This will provide a quick reference and also a place for making notations.

2. Start at the first **Test** row in Table 5-4. Under the **Scan Conditions** column highlight or type-in the CVs into the CV Name box and then enter the listed corresponding number into the Current Value box. Right mouse-click on the remaining items listed in the column.
3. Advance to the **Adjustments** column on the right and perform the directions listed in the column. Refer to the following steps for measuring the RF output power:

5-4 Head Power Monitor Functional Check Process (Continued)

- a. **If using the RF Power Measurement Kit:** Measure the amplitude of the RF waveform from the scope face in Volts peak (Vp).
- b. **If using the wattmeter procedure:** Read the wattmeter display and use the formula below to calculate the head RF power. Note that the cable loss factor was determined from the procedure in **APPENDIX D — DUMMY LOAD AND CABLES CALIBRATION**. The value calculated from the formula below is expressed in Watts and should only be compared to other values expressed in Watts.

RF Power Measurement (in watts) Using Wattmeter And Formula:
$\text{Wattmeter reading (in watts)} \times \text{cable loss factor}$

- c. **If using the oscilloscope procedure (NOT the RF Power Measurement Kit):** Read the peak voltage (Vpeak) from the scope display and use the formula below or the [Power Calculator](E:\rf\power\pwrcalc.htm) Tool (located at E:\rf\power\pwrcalc.htm on the Service Methods CD-ROM) to calculate the head RF power. Note that the dummy load and cable loss factor was determined from the procedure in **APPENDIX D — DUMMY LOAD AND CABLES CALIBRATION**. If using the Power Calculator Tool, note that the scope correction factor is equal to 1.0. The value calculated from the formula below is expressed in Watts and should only be compared to other values expressed in Watts.

RF Power Measurement (in watts) Using Oscilloscope And Formula:
$\frac{(V_{\text{peak}})^2}{100} \times \text{dummy load and cable loss factor}$

4. Confirm that the measured Vp or power meets the specifications provided in the **Verify** or **Adjustments** columns and record the PASS/FAIL status of the test in the **Verify** column in Table 5-4.
5. Perform the steps in the **Conclusion** column and then proceed to the next row, Head Peak Power Low.
6. Make the prescribed changes described in the **Scan Conditions** column in the Head Peak Power Low row.
7. Confirm that the measured Vp or power meets the specifications provided in the **Verify** or **Adjustments** columns and record the PASS/FAIL status of the test in the **Verify** column in the Head Peak Power Low row in Table 5-4.
8. Perform the steps in the **Conclusion** column and then proceed to the next row, Head Pulse Width High (Min. Limit).
9. Change the Volts/Div setting on the oscilloscope channel being used to 20 mV/Div.

5-4 Head Power Monitor Functional Check Process (Continued)

10. Make the prescribed changes described in the **Scan Conditions** column in the Head Pulse Width High (Min. Limit) row.
11. Confirm that the measured Vp or power meets the specifications provided in the **Verify** or **Adjustments** columns and record the PASS/FAIL status of the test in the **Verify** column in the Head Pulse Width High (Min. Limit) row in Table 5-4.
12. Perform the steps in the **Conclusion** column and then proceed to the next row, Head Pulse Width High (Max. Limit).
13. Make the prescribed changes described in the **Scan Conditions** column in the Head Pulse Width High (Max. Limit) row.
14. Confirm that the conditions described in the **Verify** column are satisfied and record the PASS/FAIL status of the test in the **Verify** column in the Head Pulse Width High (Max. Limit) row in Table 5-4. Note that there are no adjustments to make in the **Adjustments** column.
15. Perform the steps in the **Conclusion** column and then proceed to the next row.
16. Continue performing each of the remaining functional check rows in Table 5-4 until all are completed.

5-4 Head Power Monitor Functional Check Process (Continued)

TABLE 5-4
HEAD CHECKS (1.0T) LX

Test	Scan Conditions	Adjustments	Verify	Conclusion
HEAD PEAK POWER <i>High</i>	[Research Operations] [Display CVs] calmode = 2 trig = 7 aset = 50 [Accept] [Research Operations] [Download] [Manual Prescan]	Increase Transmit Gain (TG) until laptop or SSM LED indicates a fault occurred on one of the monitors; measure and record power or Vp as Fault 1 below. Fault 1 _____ Continue increasing TG until laptop or SSM LED indicates the other monitor has faulted; measure and record power or Vp as Fault 2 below. Fault 2 _____	Power or Vp (Fault 1, Fault 2) is within specifications: Min: _____ Vp or 160 W (52 dBm) Nom: _____ Vp or 200 W (53 dBm) Max: _____ Vp or 240 W (53.8 dBm) ☐ Pass ☐ Fail	Set TG to 0. [Done]. Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.
HEAD PEAK POWER <i>LOW</i>	[Research Operations] [Display CVs] aset = 30 [Accept] [Research Operations] [Download] [Manual Prescan]	Increase Transmit Gain (TG) until laptop or SSM LED indicates a fault occurred on one of the monitors; measure and record power or Vp as Fault 1 below. Fault 1 _____ Continue increasing TG until laptop or SSM LED indicates the other monitor has faulted; measure and record power or Vp as Fault 2 below. Fault 2 _____	Power or Vp (Fault 1, Fault 2) is within specifications: Min: _____ Vp or 100 W (50 dBm) Nom: _____ Vp or 117 W (50.7 dBm) Max: _____ Vp or 136 W (51.34 dBm) ☐ Pass ☐ Fail	Set TG to 0. [Done]. Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.

5-4 Head Power Monitor Functional Check Process (Continued)

TABLE 5-4
HEAD CHECKS (1.0T) LX (CONTINUED)

Test	Scan Conditions	Adjustments	Verify	Conclusion
HEAD PULSE WIDTH (PW) <i>High</i> (Min Limit)	[Research Operations] [Display CVs] aset = 255 calmode = 1 pwset = 35 p1 = 3912 [Accept] [Research Operations] [Download] [Manual Prescan]	Increase TG (from 0) until Vp or power measures between minimum _____ Vp or 50W (46.9 dBm) and maximum _____ Vp or 70W (48.45 dBm).	• SENSE LEDs on SSM are ON • HEAD LED on RFI is ON • FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done].
HEAD PW <i>High</i> (Max Limit)	[Research Operations] [Display CVs] p1 = 4324 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• FAULT LEDs (both ON) (Note 1) (Note 2) ☐ Pass ☐ Fail	[Done], Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.
HEAD PW <i>LOW</i> (Min Limit)	[Research Operations] [Display CVs] p1 = 1076 pwset = 10 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• SENSE LEDs on SSM are ON • HEAD LED on RFI is ON • FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done]
HEAD PW <i>LOW</i> (Max Limit)	[Research Operations] [Display CVs] p1 = 1276 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• FAULT LEDs (both ON) (Note 1) ☐ Pass ☐ Fail	[Done], Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.

Note 1

Allow at least 12 seconds after looping starts before checking status of fault LEDs.

Note 2 (Applies to Head PW High [Max Limit] ONLY)

If Fault LEDs do not illuminate, slowly increase TG in 1 unit increments until both Fault LEDs do illuminate. Verify that, at this new TG, the power level does not exceed the maximum power level previously specified in the Head Pulse Width High (Min Limit) row under the Adjustments column.

5-4 Head Power Monitor Functional Check Process (Continued)

TABLE 5-4
HEAD CHECKS (1.0T) LX (CONTINUED)

Test	Scan Conditions	Adjustments	Verify	Conclusion
HEAD DUTY CYCLE (DC) <i>High</i> (Min Limit)	[Research Operations] [Display CVs] t3 = 33333 TR_SLOP = 0 calmode = 3 pwset = 255 dcset = 130 p3 = 3900 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	- SENSE LEDs on SSM are ON - HEAD LED on RFI is ON - FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done].
HEAD DC <i>High</i> (Max limit)	[Research Operations] [Display CVs] p3 = 4767 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	FAULT LEDs (both ON) (Note 1) ☐ Pass ☐ Fail	[Done], Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.
HEAD DC <i>LOW</i> (Min Limit)	[Research Operations] [Display CVs] p3 = 750 dcset = 25 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	- SENSE LEDs on SSM are ON - HEAD LED on RFI is ON - FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done].
HEAD DC <i>LOW</i> (Max Limit)	[Research Operations] [Display CVs] p3 = 917 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	FAULT LEDs (both ON) (Note 1) ☐ Pass ☐ Fail	[Done], Place both Monitor A and B switches to RESET , then NORMAL . On laptop, press 'C' to continue.

6 - BODY HARDWARE SETUP AND CHECKS

Tests in this section allow the service engineer to selectively force body mode peak power, pulse width, and duty cycle beyond allowable limits.

6-1 Body Hardware Setup

Tests in this section allow the service engineer to selectively force body mode peak power, pulse width, and duty cycle beyond allowable limits.



PROPERTY DAMAGE! PREVENT COIL AND ASSOCIATED SWITCH DAMAGE, BY REMOVING ALL PHANTOMS AND HARDWARE (I.E., HEAD COIL, SURFACE COIL...) FROM THE MAGNET BORE.

1. Confirm that all the steps in **Section 4 - Initial Hardware Preparation** - have been done.
2. Verify that the system is not scanning and that all coils have been removed from the magnet bore. See the two **DANGER** messages on this page.



PERSONAL INJURY! PREVENT POSSIBLE RF BURNS WHEN DISCONNECTING HELIAX CABLES FROM J3 OR J4 ON THE RFI BY VERIFYING THAT THE SYSTEM IS NOT MANUALLY PRESCANNING OR SCANNING. VERIFY THAT THE SCAN DESKTOP ICON DISPLAYS THE "IDLE" MESSAGE.



3. **If using the RF Power Measurement Kit** then refer to the RF Power Measurement Kit laminated card set.
 - a. Look in the upper, right corner of each card and find the card labeled **68.8** (1.0T Body Output).



The body RF output connection is no longer to the non-existent EFB unit, as the older reference cards from the kit show, but instead to the J4 output on the rear of the RFI. An RF adapter is provided in the RF Power Measurement Kit to connect between the HN J4 body RF output and the RF Power Measurement Kit 40 dB N-connector coupler.

- b. Configure the system as in the illustration on the card.
- c. Confirm that the rotary attenuator is set to the correct position indicated on the card.
- d. Set the oscilloscope Volts/Div setting of the channel being used to 50 mV/Div.

6-1 Body Hardware Setup (Continued)

4. **If using the wattmeter or scope** (NOT the RF Power Measurement Kit): Refer to Appendix B (Alternate Equipment Setup) for the proper system body configuration.

6-2 Converting Units If Using The RF Power Measurement Kit ONLY

Perform the conversions described in this section only if using the RF Power Measurement Kit. The measurements read from the scope face will be in Volts peak (Vp). The voltage read from an accurate oscilloscope will vary from kit to kit due to variances in cable, attenuator, and directional coupler losses. It is probably easier for most people to determine the expected Volts peak-peak values based on the unique losses of the kit at this point than to wait and convert them later while the values are being read from the scope. The more experienced engineer with a good understanding of RF power measurement principles may, however, prefer to convert and determine the values as each is read from the scope face. This engineer can skip this section and instead refer to **APPENDIX C — MISCELLANIOUS RF CONVERSIONS AND FORMULAS** for the various mathematical formulas needed to perform the conversions.

1. Refer to the RF Power Measurement Kit reference card **68.8** calibration sticker (1.0T Body RF Output).
2. Record the absolute value of the **Attenuation** number (this means record it without the negative sign), listed on the calibration sticker affixed to the front of the reference card, in the specified location on Table 6-1. This is the total attenuation between the scope channel input and the J4 Body RFI RF port that is provided by the RF Power Measurement Kit. Illustration 6-1 below provides an example of what the calibration sticker looks like.



The Attenuation data in the example will be different from what is printed on the calibration sticker affixed to the card in your kit. Do not use the Attenuation data listed in the example.

1.0T BODY CARD # 68.8 CALIBRATION STICKER DATA
Card Number 46-317724P33
Nominal: 68.80 dBm
Attenuation: -65.41 dB
Scope Power: 3.39 dBm
Scope Rdg: 0.934 Vpp

SAMPLE CALIBRATION STICKER
ILLUSTRATION 6-1

6-2 Converting Units If Using The RF Power Measurement Kit ONLY (Continued)

TABLE 6-1
POWER MEASUREMENT KIT BODY ATTENUATION (IN dB)

ATTENUATION	_____dB
--------------------	---------

3. Start the LX system calculator (**NOT the RF Calculator Tool**) by moving the mouse cursor into the background, right mouse-click, and then select the calculator from the top of the Root menu list.
4. Use the LX system calculator to subtract the Attenuation value in Table 6-1 from the number in each of the 8 locations in Table 6-2.
 - a. Record the result in each of the 8 locations in Table 6-2
 - b. **Be sure to note the correct polarity sign with each answer. This is critical.** See the example column in Table 6-2.


TABLE 6-2
BODY FUNCTIONAL CHECK LIMITS (IN dBm)

	Example	Min.	Nom.	Max.
Body Peak Power High Limits in dBm	66.127 - 65.41 = <u>0.717</u>	66.127 – Attenuation = _____	66.33 – Attenuation = _____	66.53 – Attenuation = _____
Body Peak Power Low Limits in dBm	60.41 – 65.41 = <u>- 5.00</u>	60.41 – Attenuation = _____	60.79 – Attenuation = _____	61.14 – Attenuation = _____
Body Pulse Width High Limits in dBm	56.99 – 65.41 = <u>- 8.42</u>	56.99 – Attenuation = _____	N/A	58.45 – Attenuation = _____

5. Convert each of the 3 Body Peak Power High Limit values (**Min.** [Minimum], **Nom.** [Nominal], and **Max.** [Maximum]) expressed in dBm in Table 6-2 to Volts peak (Vp) using one of the three methods listed below and then record the results in the locations provided in the **Verify** column located in the first row in Table 6-4.
 - a. **METHOD 1:** Use the RF Power Calculation Tool provided on the 8.X Service Methods CD-ROM to convert from dBm to Vp.
 - i. Confirm that the 8.X Service Methods CD-ROM is in the laptop CD-ROM drive and then select the [Power Calculator](#) hyper-link to automatically start the tool.


6-2 Covertng Units If Using The RF Power Measurement Kit ONLY (Continued)

- ii. Alternatively, the tool can be started manually from Windows 95/98 by first moving the cursor to the lower left part of the screen and then selecting "**Start**" .

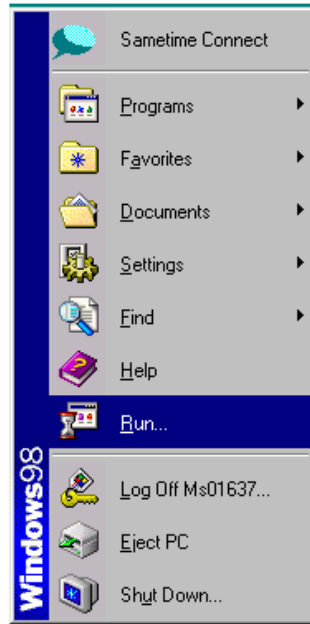

- iii. Next, select "**Run**"

- iv. Finally, type in the file pathname of **E:\rf\power\pwrcalc.htm** as shown:

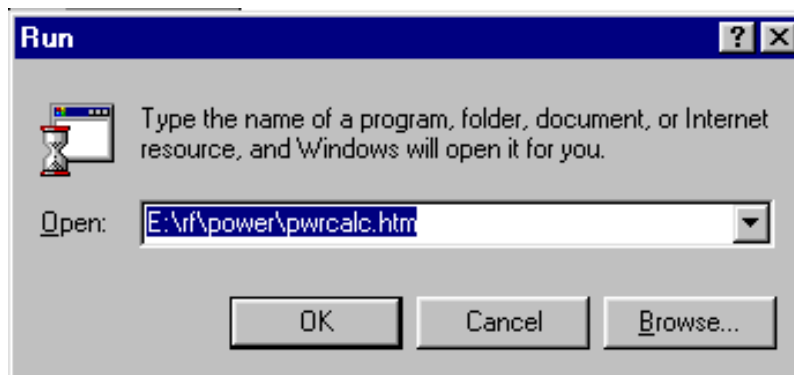

- v. Follow the directions provided in the Tool to perform the dBm to Vp conversion.
- b. **METHOD 2:** Use the LX system calculator to convert from dBm to Vp.
- i. Access the LX system calculator by moving the mouse cursor into the background, right mouse click, and then select the calculator from the Root menu list.

6-2 Coverting Units If Using The RF Power Measurement Kit ONLY (Continued)

- ii. Enter the data into the LX system calculator exactly in the order as shown below:

dBm [÷] 20 [=] [INV LOG] [*] 0.316 [=] Vp (←This is the answer)

Example: Using 0.717 dBm, the result of the first example calculation in Table 6-2, and the LX system calculator to determine the equivalent Vp:

0.717 [÷] 20 [=] 0.0359 [INV LOG] 1.086 [*] 0.316 [=] **0.343** Vp.

- c. **METHOD 3:** Use a hand calculator or the scientific calculator available on the field laptop and either of the 2 formulae listed below to convert from dBm to Vp:

Note

“antilog” may also be noted on the laptop calculator as “inv log”.

$$Vp = 10^{\left(\frac{dBm}{20}\right)} \times 0.316 \quad \text{or} \quad Vp = \text{antilog} \left(\frac{dBm}{20} \right) \times 0.316$$

where “X” signifies multiplication

6. Repeat this conversion process for the 3 Body Peak Power Low Limit values in Table 6-2 and record the results in the locations provided in the **Verify** column located in the second row in Table 6-4.
7. Repeat this conversion process for the 2 Body Pulse Width High Limit values in Table 6-2 and record the results in the locations provided in the **Adjustments** column located in the third row in Table 6-4.

6-3 Body Scan Setup

1. At operator work space, select the scan desktop ICON in the desktop control panel, if you have not already done so.
2. At the operator work space, prepare the system for a Body Power Monitor scan using the "Service Protocols" procedure in Table 6-3 below. Refer to Appendix B for the non-proprietary protocol.

TABLE 6-3
SCAN PRESCRIPTION - BODY PM CHECKS

Note: This is the alternate proprietary procedure available for GE use, and to sites with a valid Advanced Service Package Limited License. Refer to Appendix B for the non-proprietary protocol.

- A. **[New Series]**
At Patient Protocols – select **other**.
- B. In the protocol field, type **o.23.2<ENTER>** (o=Other, 23.2 =series) to load the body protocol
OR select **[o.23] [Series 2] [Accept]**.
[OK] (if required).
- C. **[Save Series]**.
- D. **[Prepare to Scan]**.
- E. **[Research Operations]**.
[Setup Params]. Set **TG** to **50**. **[Done]**.
- F. **[Research Operations]**.
[Display CVs]. Highlight CV Name and enter the following:
CV name: **dcset <ENTER>, 255 <ENTER>**
CV name: **t3 <ENTER>, 20000 <ENTER>**
[Accept].

6-4 Body Power Monitor Functional Check Process

1. Refer to Table 6-4. Note that the table consists of 5 columns (**Test, Scan Conditions, Adjustments, Verify, Conclusion**) and 10 individual rows of various head power tests.

Note

It may be advantageous at this point to print a hard copy of Table 6-4. This will provide a quick reference and also a place for making notations.

2. Start at the first **Test** row in Table 6-4. Under the **Scan Conditions** column highlight or type-in the CVs into the CV Name box and then enter the listed corresponding number into the Current Value box. Right mouse-click on the remaining items listed in the column.
3. Advance to the **Adjustments** column on the right and perform the directions listed in the column. Refer to the following steps for measuring the RF output power:

6-4 Body Power Monitor Functional Check Process (Continued)

- a. **If using the RF Power Measurement Kit:** Measure the amplitude of the RF waveform from the scope face in Volts peak (Vp).
- b. **If using the wattmeter procedure:** Read the wattmeter display and use the formula below to calculate the body RF power. Note that the dummy load and cable loss factor was determined from the procedure in **APPENDIX D — DUMMY LOAD AND CABLES CALIBRATION**. The value calculated from the formula below is expressed in Watts and should only be compared to other values expressed in Watts.

RF Power Measurement (in watts) Using Wattmeter And Formula:

Wattmeter reading (in watts) X dummy load and cable loss factor
--

- c. **If using the oscilloscope procedure (NOT the RF Power Measurement Kit):** Read the peak voltage (Vpeak) from the scope display and use the formula below or the [Power Calculator](E:\rf\power\pwrcalc.htm) Tool (located at E:\rf\power\pwrcalc.htm on the Service Methods CD-ROM) to calculate the body RF power. Note that the dummy load and cable loss factor was determined from the procedure in **APPENDIX D — DUMMY LOAD AND CABLES CALIBRATION**. If using the Power Calculator Tool, note that the scope correction factor is equal to 1.0. The value calculated from the formula below is expressed in Watts and should only be compared to other values expressed in Watts.

RF Power Measurement (in watts) Using Oscilloscope And Formula:
--

$\frac{(V_{peak})^2}{100} \times \text{dummy load and cable loss factor}$

4. Confirm that the measured Vp or power meets the specifications provided in the **Verify** or **Adjustments** columns and record the PASS/FAIL status of the test in the **Verify** column in Table 6-4.
5. Perform the steps in the **Conclusion** column and then proceed to the next row, Body Peak Power Low.
6. Make the prescribed changes described in the **Scan Conditions** column in the Body Peak Power Low row.
7. Confirm that the measured Vp or power meets the specifications provided in the **Verify** or **Adjustments** columns and record the PASS/FAIL status of the test in the **Verify** column in the Body Peak Power Low row in Table 6-4.
8. Perform the steps in the **Conclusion** column and then proceed to the next row, Body Pulse Width High (Min. Limit).
9. Change the Volts/Div setting on the oscilloscope channel being used to 20 mV/Div.

6-4 Body Power Monitor Functional Check Process (Continued)

10. Make the prescribed changes described in the **Scan Conditions** column in the Body Pulse Width High (Min. Limit) row.
11. Confirm that the measured Vp or power meets the specifications provided in the **Verify** or **Adjustments** columns and record the PASS/FAIL status of the test in the **Verify** column in the Body Pulse Width High (Min. Limit) row in Table 6-4.
12. Perform the steps in the **Conclusion** column and then proceed to the next row, Body Pulse Width High (Max. Limit).
13. Make the prescribed changes described in the **Scan Conditions** column in the Body Pulse Width High (Max. Limit) row.
14. Confirm that the conditions described in the **Verify** column are satisfied and record the PASS/FAIL status of the test in the **Verify** column in the Body Pulse Width High (Max. Limit) row in Table 6-4. Note that there are no adjustments to make in the **Adjustments** column.
15. Perform the steps in the **Conclusion** column and then proceed to the next row.
16. Continue performing each of the remaining functional check rows in Table 6-4 until all are completed.

6-4 Body Power Monitor Functional Check Process (Continued)

TABLE 6-4
BODY CHECKS (1.0T) LX

Test	Scan Conditions	Adjustments	Verify	Conclusion
BODY PEAK POWER <i>HIGH</i>	[Research Operations] [Display CVs] calmode = 2 trig = 7 aset = 110 [Accept] [Research Operations] [Download] [Manual Prescan]	Increase Transmit Gain (TG) until laptop or SSM LED indicates a fault occurred on one of the monitors; measure and record power or Vp as Fault 1 below. Fault 1 _____ Continue increasing TG until laptop or SSM LED indicates the other monitor has faulted; measure and record power or Vp as Fault 2 below. Fault 2 _____	Power or Vp (Fault 1, Fault 2) is within specifications: Min: _____ Vp or 4100W (66.127dBm) Nom: _____ Vp or 4300W (66.33dBm) Max: _____ Vp or 4500W (66.53dBm) ☐ Pass ☐ Fail	Set TG to 0. [Done] Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.
BODY PEAK POWER <i>LOW</i>	[Research Operations] [Display CVs] aset = 30 [Accept] [Research Operations] [Download] [Manual Prescan]	Increase Transmit Gain (TG) until laptop or SSM LED indicates a fault occurred on one of the monitors; measure and record power or Vp as Fault 1 below. Fault 1 _____ Continue increasing TG until laptop or SSM LED indicates the other monitor has faulted; measure and record power or Vp as Fault 2 below. Fault 2 _____	Power or Vp (Fault 1, Fault 2) is within specifications: Min: _____ Vp or 1100W (60.41dBm) Nom: _____ Vp or 1200W (60.79dBm) Max: _____ Vp or 1300W (61.14dBm) ☐ Pass ☐ Fail	Set TG = 0. [Done] Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.

6-4 Body Power Monitor Functional Check Process (Continued)

TABLE 6-4
BODY CHECKS (1.0T) LX (CONTINUED)

Test	Scan Conditions	Adjustments	Verify	Conclusion
BODY PULSE WIDTH (PW) HIGH (Min Limit)	[Research Operations] [Display CVs] aset = 255 calmode = 1 pwset = 35 p1 = 3912 [Accept] [Research Operations] [Download] [Manual Prescan]	Increase TG (from 0) until Vp or power measures between minimum _____ Vp or 500W (56.99 dBm) and maximum _____ Vp or 700W (58.45 dBm).	• SENSE LEDs on SSM are ON, • BODY LED on RFI is ON, • FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done]
BODY PW HIGH (Max Limit)	[Research Operations] [Display CVs] p1 = 4324 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• FAULT LEDs (both ON) (Note 1) (Note 2) ☐ Pass ☐ Fail	[Done] Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.
BODY PW LOW (Min Limit)	[Research Operations] [Display CVs] p1 = 1076 pwset = 10 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• SENSE LEDs on SSM are ON, • BODY LED on RFI is ON, • FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done]
BODY PW LOW (Max Limit)	[Research Operations] [Display CVs] p1 = 1276 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• FAULT LEDs (both ON) (Note 1) ☐ Pass ☐ Fail	[Done] Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.

Note 1

Allow at least 12 seconds after looping starts before checking status of fault LEDs.

Note 2 (Applies to Body PW High [Max Limit] ONLY)

If Fault LEDs do not illuminate, slowly increase TG in 1 unit increments until both Fault LEDs do illuminate. Verify that, at this new TG, the power level does not exceed the maximum power level previously specified in the Body Pulse Width High (Min Limit) row under the Adjustments column.

6-4 Body Power Monitor Functional Check Process (Continued)

TABLE 6-4
BODY CHECKS (1.0T) LX (CONTINUED)

Test	Scan Conditions	Adjustments	Verify	Conclusion
BODY DUTY CYCLE (DC) <i>HIGH</i> (Min Limit)	[Research Operations] [Display CVs] calmode = 3 t3 = 33333 TR_SLOP = 0 pwset = 255 dcset = 130 p3 = 3900 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• SENSE LEDs on SSM are ON, • BODY LED on RFI is ON, • FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done]
BODY DC <i>HIGH</i> (Max Limit)	[Research Operations] [Display CVs] p3 = 4767 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• FAULT LEDs (both ON) (Note 1) ☐ Pass ☐ Fail	[Done] Place both Monitor A and B switches to RESET , then BYPASS . On laptop, press 'C' to continue.
BODY DC <i>LOW</i> (Min Limit)	[Research Operations] [Display CVs] p3 = 750 dcset = 25 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• SENSE LEDs on SSM are ON, • BODY LED on RFI is ON, • FAULT LEDs (both OFF) (Note 1) ☐ Pass ☐ Fail	[Done]
BODY DC <i>LOW</i> (Max Limit)	[Research Operations] [Display CVs] p3 = 917 [Accept] [Research Operations] [Download] [Manual Prescan]	Do not change TG.	• FAULT LEDs (both ON) (Note 1) ☐ Pass ☐ Fail	[Done] Place both Monitor A and B switches to RESET , then NORMAL . On laptop, press 'C' to continue.

7 - AMPLIFIER SHUTDOWN VERIFICATION

Tests in this section verify that each power monitor is capable of shutting down the RF Amplifier independently of the other, and further verify both the logic and relay shutdowns for each channel. Setpoints are not checked in this section — this test is forcing a shutdown.

When one or both of the monitors is in bypass mode, the HV relay shutdown is disabled.

1. Enable power monitor A by placing the monitor A switch on the front panel of the SSM in the Normal (down) position.
2. At the operator work space, prepare the system for a Amp Shutdown Verification scan using the "Service Protocols" procedure in Table 7-1 below. Refer to the body protocol in Appendix B for the non-proprietary protocol.

Note

Verify that the Head Coil has been removed from the patient table before continuing.

TABLE 7-1

SCAN PRESCRIPTION - HIGH VOLTAGE RELAY SHUTDOWN — BODY MODE

Note: This is the alternate proprietary procedure available for GE use, and to sites with a valid Advanced Service Package Limited License.

A. **[New Series]**

At Patient Protocols – select **other**.

- B. In the protocol field, type **o.23.3<ENTER>** (o=Other, 23.3 =series) to load the body protocol
OR select **[o.23] [Series 3] [Accept]**.
[OK] (if required).

C. **[Save Series]**

D. **[Prepare to Scan]**.

E. **[Research Operations]**.

[Setup Params]. Set **TG** to **0**. **[Done]**.

F. **[Research Operations]**.

[Display CVs]. Highlight CV Name and enter the following:

CV name: **trig <ENTER>, 1 <ENTER>**

CV name: **aset <ENTER>, 120 <ENTER>**

CV name: **calmode <ENTER>, 2 <ENTER>**

CV name: **pwset <ENTER>, 255 <ENTER>**

CV name: **dcset <ENTER>, 255 <ENTER>**

CV name: **p1 <ENTER>, 3100 <ENTER>**

CV name: **p3 <ENTER>, 2400 <ENTER>**

G. **[Accept]**

H. **[Research Operations]**

I. **[Download]**.

J. Select **[Manual Prescan]**.

7 - AMPLIFIER SHUTDOWN VERIFICATION (CONTINUED)



The system error log will be checked in the next step for a power monitor fault error message. Make sure the error log Viewing Level selection is set to “ALL” before reading the error log. The error message may not otherwise be visible.

3. Increase TG until monitor A FAULT LED illuminates. Quickly verify the following:
 - The FAULT LED on for monitor A (resets automatically after 1 second).
 - Software Message Log - Software will log a power monitor fault in the message log.
 - Verify oscilloscope or Wattmeter (press FWD PEP then MAX) reads zero RF power.

Note

Do not perform any operations at the keyboard while In Place Error Recovery (IPER) is resetting the RF amplifier.

4. While waiting for the system to recover, place monitor B in normal mode by placing the monitor B switch on the front panel of the SSM in the Normal (down) position. Now both monitors A and B should be in normal mode.
5. After IPER has reset the RF amplifier (the “Please Press Start Scan Button” message appears in the upper left message window), decrease TG by 25 counts (2.5 dB) and then press START SCAN on the operator's console to initiate looping.
6. Increase TG until the FAULT LED for monitor A or B illuminates. Quickly verify the following:
 - FAULT LED on for monitor A or B (resets automatically after two seconds).
 - Software Message Log - Software will log a power monitor fault in the message log.
 - Verify oscilloscope or Wattmeter (press FWD PEP then MAX) reads zero RF power.

Note

Do not perform any operations at the keyboard while In Place Error Recovery (IPER) is resetting the RF amplifier.

7. While waiting for the system to recover, disable power monitor A by placing the monitor A switch on the front panel of the SSM in the Bypass position.
8. After IPER has reset the RF amplifier (the “Please Press Start Scan Button” message appears in the upper left message window), decrease TG by 25 counts (2.5 dB) and then press START SCAN on the operator's console to initiate looping.
9. Increase TG until monitor B FAULT LED illuminates. Quickly verify the following:
 - FAULT LED on for monitor B (resets automatically after two seconds).
 - Software Message Log - Software will log a power monitor fault in the message log.
 - Verify oscilloscope or Wattmeter (press FWD PEP then MAX) reads zero RF power.

7 - AMPLIFIER SHUTDOWN VERIFICATION (CONTINUED)


Note

Do not perform any operations at keyboard while In Place Error Recovery (IPER) is resetting the RF amplifier.

10. After IPER has reset the RF amplifier (the "Please Press Start Scan Button" message appears in the upper left message window), select **[Done]**.
11. In the Rx Manager, select **[End Exam]**, then **[Confirm]**.
12. Proceed to **SECTION 8 – SYSTEM RESTORATION**.

8 - SYSTEM RESTORATION

1. Verify the system is not scanning and the scan desktop icon is displaying the "Idle" message.
2. Refer to Illustration 8-1. At the front of the SSM place the:
 - 2 (two) power MONITOR switches (A and B) to the bottom NORMAL position.
 - 3 (three) DRV FLT switches to the top EN (enable faults) position.

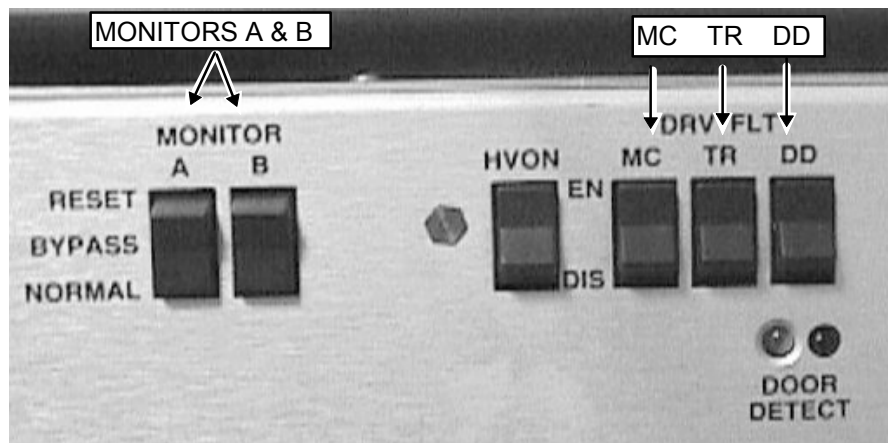

FRONT PANEL SWITCHES ENABLED
ILLUSTRATION 8-1

3. Place DIP switches 4, 3, 2, and 1 of the Test Switches in the DOWN position or disconnect laptop computer from J46 on the rear of the SSM.

Note

If using the Laptop perform the following steps to exit MONS1.EXE and return to Windows:

A) Follow screen directions to return to DOS.

B) At the DOS C:\CCLASS\RFTOOLS> prompt, type **EXIT<ENTER>**.

4. Ensure that the dummy load (all test equipment) is disconnected from the rear of the cabinet.
5. Ensure that the head and body Heliax cables are reconnected to the rear of the cabinet.
6. Complete one head scan satisfactorily.
7. Complete one body scan satisfactorily.
8. Replace the front and rear covers on the cabinet.

Prepare and file power monitor functional check data sheet according to current company policy. See **SECTION 9 – DECLARATION FORM PREPARATION**.

9 - DECLARATION FORM PREPARATION

1. Prepare Direction 2212504, *Signa Release 8.X RF/PDU 1.0T/1.5T Cabinet Power Monitor Functional Test Declaration Form, Appendix A*. Refer to procedure for Data Sheet for Power Monitor Test Declaration Form.
2. File completed declaration form per current company policy.

APPENDIX A — SIGNA POWER MONITOR FUNCTIONAL TEST DECLARATION FORM

Signa Release 8.X RF/PDU 1.0T/1.5T CABINET Power Monitor Functional Test Declaration

SITE

(NAME) _____ (SYSTEM ID, ROOM NUMBER, ETC., IF APPLICABLE) _____
 (STREET) _____ (SYSTEM CABINET SERIAL NUMBER) _____
 (CITY, STATE, ZIP CODE) _____ (CURRENT SYSTEM SOFTWARE RELEASE) _____

TEST EQUIPMENT

(BIRD WATTMETER SERIAL NO.) _____ (OSCILLOSCOPE SERIAL NO.) _____ (RF POWER MEASUREMENT KIT SERIAL NO.) _____
 (DUMMY LOAD PLUS CABLES ATTEN FACTOR) _____ (AMP TO WATTMETER CABLE ATTEN FACTOR) _____

PERFORMANCE

FUNCTION CHECK STATUS	PASS
BODY PEAK POWER TEST	☐
BODY PULSE WIDTH TEST	☐
BODY DUTY CYCLE TEST	☐
HEAD PEAK POWER TEST	☐
HEAD PULSE WIDTH TEST	☐
HEAD DUTY CYCLE TEST	☐
SHUTDOWN VERIFICATION	☐

CONFIGURATION FILE ENTRIES		
	Cable loss coefficient	Coil loss coefficient
Head Coil	_____	_____
Body Coil	_____	_____
AMP Cal Head	_____	ABSOLUTE AVG SAR _____
AMP Cal Body	_____	ABSOLUTE PEAK SAR _____
Average SAR Factor	_____	Peak SAR Factor _____

RESTORATION CHECK LIST

- ☐ On the front panel of the SSM, place the two monitor switches in the Normal position, and place the three DRV FLT switches in the enable faults position (EN).
- ☐ Ensure that the dummy load is disconnected from the rear of the cabinet.
- ☐ Ensure that the head and body Helix cables are reconnected to the rear of the cabinet.
- ☐ Disconnect laptop computer from J46 on the rear of the SSM.
- ☐ Complete one head scan satisfactorily.
- ☐ Complete one body scan satisfactorily.
- ☐ Replace the front door on the cabinet.
- ☐ Prepare and file Direction 2212504, *Signa Horizon 5.X & 8.X RF/PDU 1.0T/1.5T Power Monitor Functional Test Declaration Form* according to company policy.

(PERFORMED BY - PLEASE PRINT)

(SIGNATURE)


(DATE)

APPENDIX B — ALTERNATE EQUIPMENT SETUPS & NON-SERVICE PROTOCOLS

Use this section if the RF Power Measurement Kit is NOT going to be used. This section contains the setup information and protocols necessary to measure the body and head RF output power using either a wattmeter or oscilloscope.

TABLE B-1
TOOLS AND INSTRUMENTS REQUIRED WITHOUT RF POWER MEASUREMENT KIT

Item	Description	Part Number
1.	RF Test Cables Kit	46-255816G1
2.	50 ohm, 200 Watt, 30 dB attenuator - Bird Model 8322	46-255837P10
3.	Oscilloscope	46-183029P61
4.	Laptop serial cable (optional)	2124497-46
4.	Wattmeter and appropriate elements (optional)	Not supplied

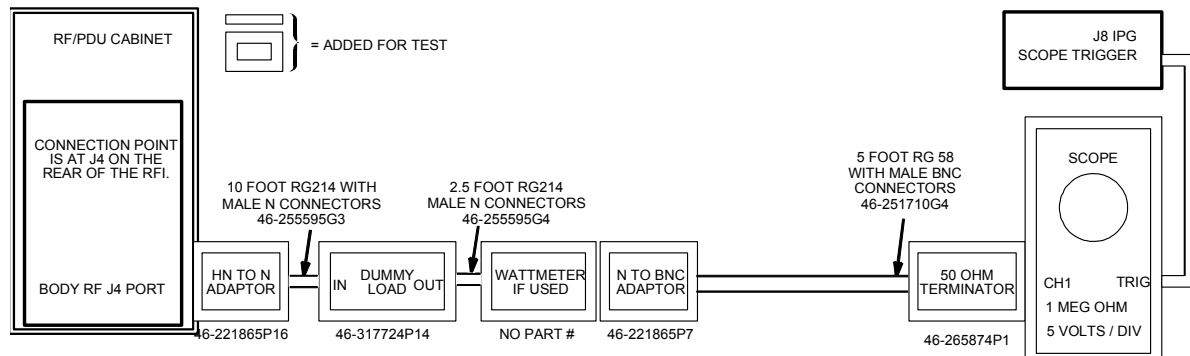

ALTERNATE BODY EQUIPMENT SETUP FOR WATTMETER OR SCOPE MEASUREMENT
ILLUSTRATION B-1

Note

If the wattmeter is not available then bypass it's connection in the circuit using an N-female to N-female connector 46-265875P2.


NON-SERVICE BODY PROTOCOL

<u>PATIENT REGISTER</u>	[New Pt]
<u>PATIENT INFORMATION</u>	
Patient Id	geservice
Patient Name	body pmc
Weight (Lb)	300 - IMPORTANT
Landmark	[Landmark] [>] [Sternal Notch]
<u>PATIENT PROTOCOLS</u>	[Patient Position]
<u>PATIENT POSITION</u>	
Patient Position	[>] [Supine]
Patient Entry	[>] [Head First]
Coil	[...] [Body] [Accept]
<u>IMAGING PARAMETERS</u>	
Plane	[>] [Axial]
Mode	[>] [2D]
Pulse Seq	[...] [Spin Echo] [Accept]
Imaging Options	none (default)
Psd Name	cal
Protocol	no entry

<u>SCAN TIMING</u>	
* of Echoes	1 (default)
TE	[25]
TR	[200]

<u>SCANNING RANGE</u>	
FOV	[24]
Slice Thickness	[10]
Spacing	0
Start	0
End	0
# Slices	1 (default)
L/R Center	0 (default)
P/A Center	0 (default)
Table Delta	0.00 (default)
<u>ACQUISITION TIMING</u>	
Freq	[256]
Phase	[128]
NEX	[2]
Freq Dir	[>] [R/L]
Auto Center Freq	[>] [Peak]
(lowest window)	[Save Series]
[Research Operations]	[Setup Params]
Set TG	50 [Done]
[Research Operations]	[Download]
[Prepare to Scan]	

If performing Amplifier Shutdown Verification of Section 7 then be sure to enter the CVs per Table 7-1, items F - H.

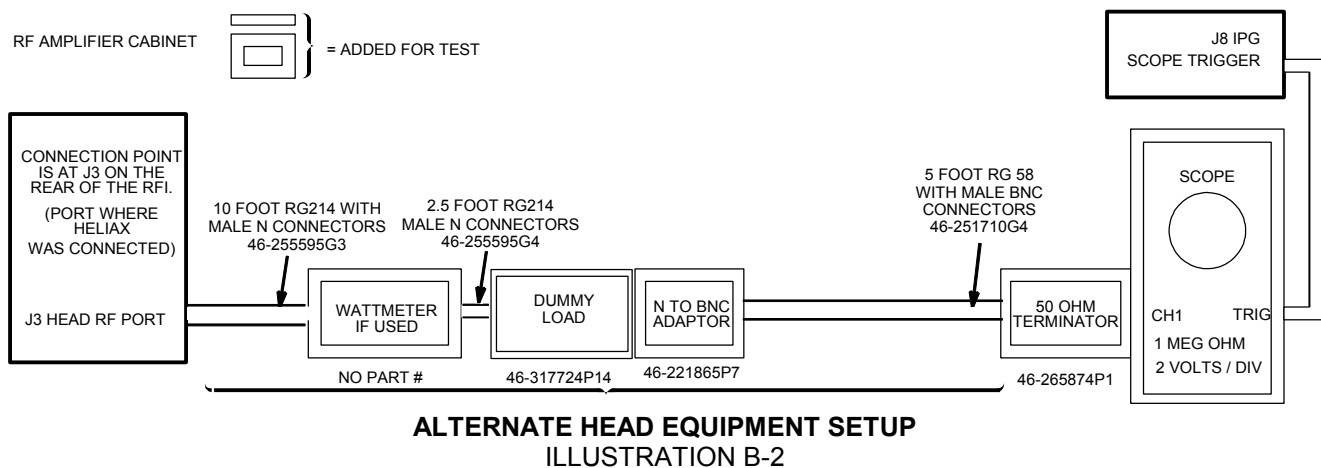

Note

If the wattmeter is not available then bypass it's connection in the circuit using an N-female to N-female connector 46-265875P2.

NON- SERVICE HEAD PROTOCOL

<u>PATIENT REGISTER</u>	[New Pt]	<u>SCANNING RANGE</u>
<u>PATIENT INFORMATION</u>		FOV [24]
Patient Id	geservice	Slice Thickness [10]
Patient Name	head pmc	Spacing 0
Weight (Lb)	300 - IMPORTANT	Start 0
	[Landmark]	End 0
Landmark	[>] [Sternal Notch]	# Slices 1 (default)
		L/R Center 0 (default)
<u>PATIENT PROTOCOLS</u>	[Patient Position]	P/A Center 0 (default)
		Table Delta 0.00 (default)
<u>PATIENT POSITION</u>		<u>ACQUISITION TIMING</u>
Patient Position	[>] [Supine]	Freq [256]
Patient Entry	[>] [Head First]	Phase [128]
Coil	[...] [Head] [Accept]	NEX [2]
<u>IMAGING PARAMETERS</u>		Freq Dir [>] [R/L]
Plane	[>] [Axial]	Auto Center Freq [>] [Peak]
Mode	[>] [2D]	(lowest window) [Save Series]
Pulse Seq	[...] [Spin Echo]	[Research Operations] [Setup Params]
	[Accept]	Set TG 50 [Done]
Imaging Options	none (default)	[Research Operations] [Download]
Psd Name	cal	[Prepare to Scan]
Protocol	no entry	
<u>SCAN TIMING</u>		
* of Echoes	1 (default)	
TE	[25]	
TR	[200]	

APPENDIX C — SOFTWARE AND FORMULAE FOR MISCELLANEOUS RF CONVERSIONS

Software are described and formulae are presented in this section that will allow the FE to perform various conversions and calculations.

C-1 RF Conversion Software

C-1-1 RF Power Calculator

The first tool, the RF Power Calculator, is a UNIX-based application available on the LX system under the **Toolbelt icon, [Utilities], [RF Calculator], [Start]**. Watts to dBm, dBm to Watts, relative volts to relative power, and relative power to relative volts conversions can be done with the RF Power Calculator tool.

C-1-2 Power Calculator Tool

The second tool, the [Power Calculator](#), is an HTML document that is available from the 8.X Service Methods CD-ROM. This tool will run from either the laptop PC or MR host computer. The tool can manually be started from the laptop PC by selecting **Start, Run**, and then entering the following file pathname in the “Open” entry box: E:\rf\power\pwrcalc.htm. Next, select **OK**. The tool is divided into 3 sections. The first section allows one to calculate RF power from measured peak voltage when using a 100 MHz oscilloscope and a properly characterized RF cables kit and dummy load. Note that at 1.0T the scope correction factor is equal to 1.0. The second section provides Vpp to dBm, dBm to Watts, and dBm to Vpp conversion capability. The third section provides short instructions and a means for comparing the attenuation in dB reported on each RF Power Measurement Kit Card with the measured Magnitude Squared Attenuation Factor reported by the Attenuation Test tool that is described in **APPENDIX D – DUMMY LOAD AND CABLES CALIBRATION**.

C-2 RF Conversion Formulae

These are provided to assist the experienced FE with RF power measurement and troubleshooting. The RF conversions below are presented in two different formats. The first format is the actual mathematical formula. The second format assumes that you are entering the data into the LX system calculator exactly as you read it (Access the LX system calculator by right mouse clicking in the background and then selecting the calculator from the Root menu.).

VP-P to dBm Calculation:

$$\text{dBm} = 20 \log \left(\frac{V_{pp}}{0.632} \right)$$

VP-P [÷] 0.632 [=] [LOG] [*] 20 [=] dBm.

APPENDIX C — SOFTWARE AND FORMULAE FOR MISCELLANEOUS RF CONVERSIONS (CONTINUED)

dBm to VP-P Calculation:

$$V_{pp} = 10^{\left(\frac{dBm}{20}\right)} \times 0.632 \quad \text{or} \quad V_{pp} = \text{antilog} \left(\frac{dBm}{20} \right) \times 0.632$$

dBm [+] 20 [=] [INV LOG] [*] 0.632 [=] VP-P.

Example: Using 3.65 dBm, the Scope Power reading listed on a particular Kit card, and the LX system calculator to determine the equivalent Vpp:

3.65 [+] 20 [=] (0.1825) [INV LOG] (1.5222991) [*] 0.632 [=] (0.9620931) VP-P.

dBm to Watts Calculation:

$$\text{Watts} = 10^{\left(\frac{dBm}{10}\right)} \times 0.001 \quad \text{or} \quad \text{Watts} = \text{antilog} \left(\frac{dBm}{10} \right) \times 0.001$$

dBm [+] 10 [=] [INV LOG] [=] [*] 0.001 [=] Total Watts.

APPENDIX D — DUMMY LOAD AND CABLES CALIBRATION

Description - This procedure provides directions for determining the true loss attributable to the dummy load and cables used when measuring the RF output power with a wattmeter or oscilloscope of 100 MHz bandwidth or greater. It is necessary to know and account for the actual loss these components contribute in order to accurately measure RF power. This procedure is not needed if using the RF Power Measurement Kit. The RF Power Measurement Kit has already been calibrated so that this loss is known and accounted for.

Note

This tool will accurately measure attenuation values up to 50dB. The combined attenuation of the dummy load and cables is well below this value. The accuracy of the tool begins to diminish when measuring attenuation values over 50dB.

D-1 Overview

Test cables long enough to reach the cables, connectors and dummy load to be tested are connected between the Exciter RF Output (System Cabinet J1) and Receiver Body Input (System Cabinet J2). Receiver gains (R1 & R2) and transmit gain (TG) are set for near full scale reading on the power spectrum during prescan calibration. A reference scan is taken and stored in a raw file. The Attenuation Test Tool is used to calculate the baseline factor from the reference scan for the test cables (i.e., there is some loss from the test cables).

The dummy load and cable(s) to be tested are next inserted in series with the test cables and another scan is taken. Again, the Attenuation Test Tool is used to determine the “Magnitude Squared Attenuation Factor” (i.e., how much has the test signal been attenuated?). This attenuation factor is used in the RF power calibration process to accurately calculate the RF power level.

Note

If any problems are encountered during the following procedure, always start over at the beginning and re-do the reference scan. Then you may add, as directed in this procedure, any type of attenuation hardware you might have reason to test.

D-2 Tools and Instruments Required

See Table D2-1

TABLE D2-1
REQUIRED TOOLS AND INSTRUMENTS

Item	Description	Part Number	Qty.
1	50-ohm dummy load, 200 watt, 30 dB attenuator - Bird Model 8322 (or equivalent).	46-317724P14	1
2	RF Test Cables Kit	46-255816G1	1


D-3 Initial Setup

1. Open the System Cabinet back door and locate the TNS module inside the upper left of the cabinet. This unit has an LCD display, Reset Button, and Disable/Enable toggle switch on the front.
2. Locate the body TNS on the TNS module. This is one of the two shiny metal TNS boxes affixed to the main TNS module assembly that is mounted farthest from the rear of the cabinet. Multicoil systems have an additional piggyback board affixed with 3 extra TNS boxes that mounts over the top of the main TNS module assembly.
3. Bypass the body TNS (A24 A1 A1) from the circuit by disconnecting the small coax cables A24 A1 A1 J12 (signal output) and A24 A1 A1 J2 (signal input) from the TNS and connecting both together using a female-BNC to female-BNC adapter (also known as a bullet adapter 46-220427P3).

Note

Failure to bypass the body TNS out of the circuit may result no signal being received by the system. Merely disabling the TNS by moving the Enable/Disable toggle switch down to the Disable position instead of bypassing it out of the circuit may work, however, the TNS processor can override the switch.

4. Reconfigure test hardware as shown in Illustration D3-1.

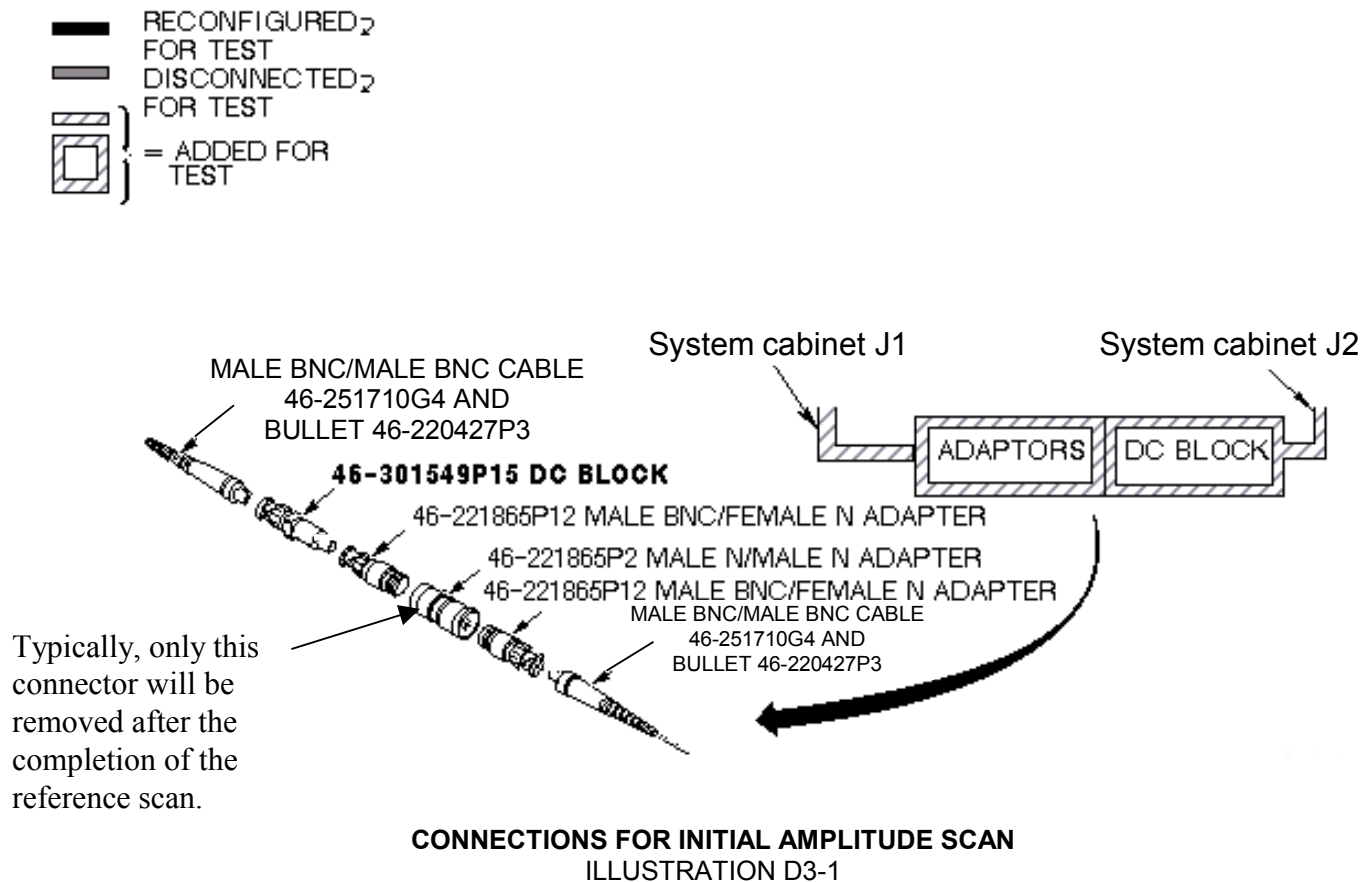

Note

Adapters not shown in Illustration D3-1 can be added, if necessary, from the RF Cables Kit. Usage of the inline DC Block as shown in Illustration D3-1 is mandatory.

5. Disconnect the existing cables at the Systems cabinet I/F panel J1 (RF Out) and J2 (Receiver Body Input) and set them aside.
6. Connect the assembled test cables and adapters between J1 (RF Out) and J2 (Receiver Body Input) on the System Cabinet interface.

7. At the operator workspace, prepare the system for a Dummy Load scan using the procedure, see below.
 - a. Click on **[New Pt]**
 Id: **geservice**
 Name: **dummy load**
 Weight (Lb): **111**
 Set Patient Protocols to **Service**.
 - b. In the Protocol field, type **o.18.1** (o=Other, 1=series number) to load the protocol.
 - c. Set a landmark if necessary, then **Save Series**.
 - d. With the right mouse key, select **[Research Operations]**, then select **[Display CVs]**.
 Set value of CV **calmode** to **2** (trapezoid pulse).
 (Caution here. Make sure the previous CV has been cleared before entering the next one. Look at the screen!)
 Set value of CV **p2_ramp** to **1** (1 μ sec ramp time).
 Set value of CV **t2** to **50000** (50 msec tr).
 Set value of CV **pismode** to **1** (exc service).
 Set value of CV **pmode** to **1** (data collection).
 Set value of CV **daqm** to **1** (data in window).
 - e. Select **[Accept]** and then select **[Research Operations]**, then select **[Download]** then select **[Manual Prescan]**.

D-4 Data Collection

1. When in **[Manual Prescan]**, set **R1** to **7**, and **R2** to **14**.
2. Adjust transmit gain (**TG**) to achieve an R1 or R2 (on IP display) of approximately 98%, without going over.
3. Select **[Done]**.
4. Select **[Scan]** (Ignore the message: MR signal too large, reduce receiver gain.) (Note: on the LX systems tested, the scan time starts at 13 seconds, counts down to 7 seconds, then ends. This is normal and is not cause for alarm.)
5. From the MR Tools desktop, select **[Cals/checks]** and then **[Attenuation Test]**.
6. Use **[Atten Test]** tool selection to analyze data, as shown in Table D4-1.

TABLE D4-1
DATA COLLECTION

Output/Prompts	Input/Comments
Last run number used was: XXXX Please enter runfile number (XXXX):	<Enter>
Please select Locked / Unlocked file (L,U) (U):.....	<Enter> (working)
<pre> ***** ***** ***** Average Max. magnitude Across All Views = aaaaa Average Max. magnitude Squared = bbbbb Average RMS Across All Views = ccccc ***** ***** ***** </pre>	
Do you want to make this run the reference(Y,N)(N):.....	Y<Enter>
<i>STOP!</i> Do not answer the next question at this point. Continue with step 7 below.	


7. The next step will involve removing only the center male N to male N adapter from the test cables, setting it aside, and adding in the dummy load and cables that need to be characterized.

Note


This assumes that the item to be characterized has N female connectors at it's input and output. If it has an N connector at the input and a BNC at the output or BNC connectors at both the input and output then an additional adapter(s) will be needed in order to connect it to the test cables. Adding in one or two uncharacterized adapters should not appreciably change the baseline attenuation factor. In this case, if the wattmeter is not being used, the 2.5 ft. cable (RG214, 46-255595G4) can be eliminated from the circuit. It is not needed.

8. Connect your test cables to the opposite ends of either:

- Illustration D4-1 - Dummy Load and Cables
- Illustration D4-2 - Amplifier to Wattmeter Cable.

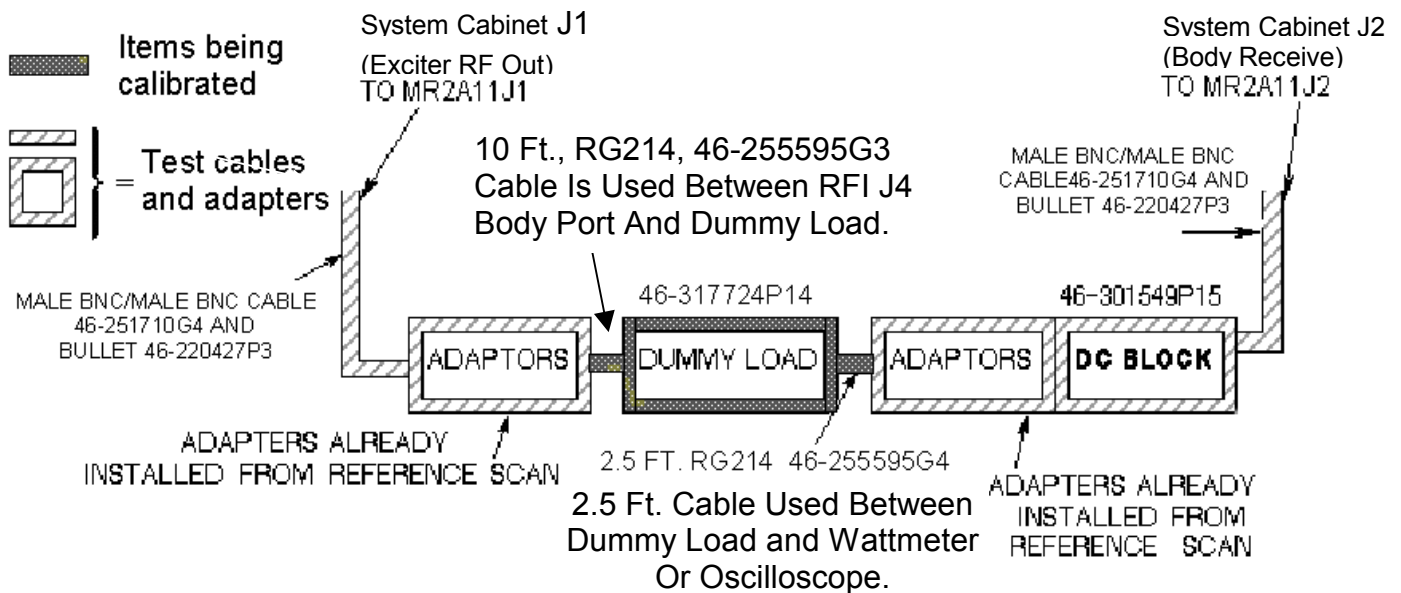

CONNECTIONS FOR DUMMY LOAD + CABLES SCAN
ILLUSTRATION D4-1

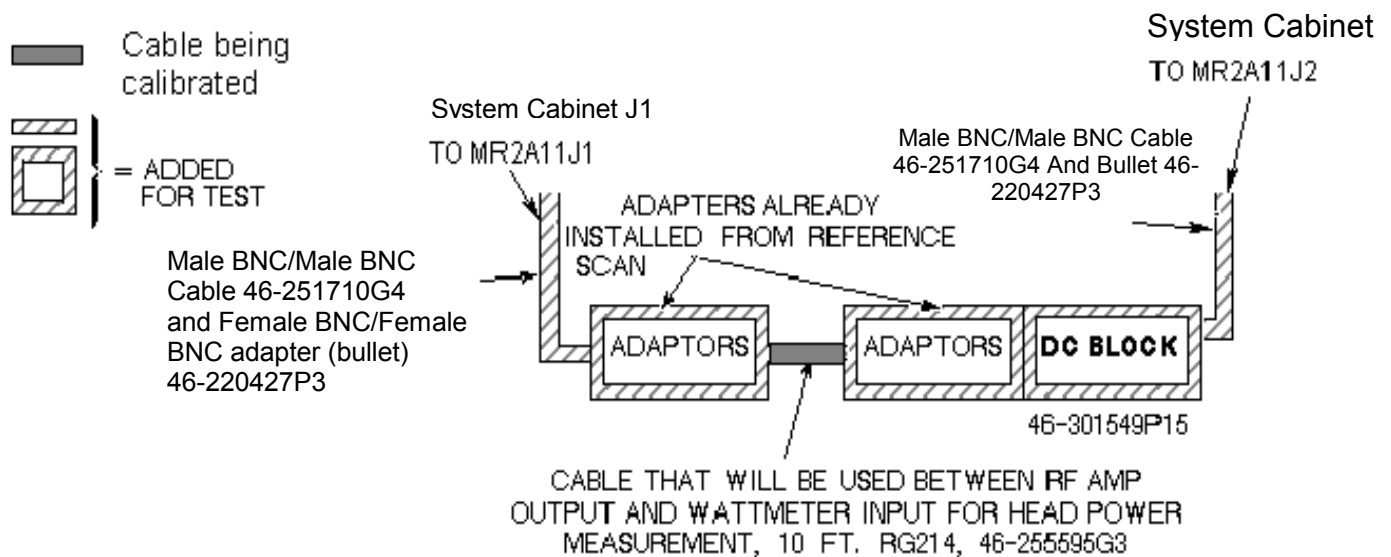

AMP TO WATTMETER CABLE SCAN
ILLUSTRATION D4-2

Note

"Bullet" RF connector referred to in illustrations D4-1 and D4-2 is a female BNC/female BNC adapter, 46-220427P3.

9. Select the scanning icon again to activate the scanning screen.
10. Select **[Scan]**.
11. When the scan is completed, re-select the tools icon again and begin the Analysis in section D-5. It will be necessary to toggle between the Scan and Toolbelt icons if multiple passes are done.

D-5 Analysis

See Table D5-1.

TABLE D5-1
ANALYSIS

Output/Prompts	Input/Comments
Do you want to compute Gain or Attenuation Ratio(G,A)[G]:	A<Enter> IMPORTANT: Answer after the scan is done.
Last run number used was: XXXX	
Please enter runfile number	<Enter>
[XXXX]:.....	
Please select Locked / Unlocked file (L,U)	<Enter> (Working)
[U]:.....	

Average Max. Magnitude Across All Views = aaaaa	
Average Max. Magnitude Squared = bbbbb	
Average RMS Across All Views = ccccc	
Magnitude Attenuation Factor = xxxxx	
Magnitude Squared Attenuation Factor = yyyyy	<=====Record in "Value" column in Table D5-2.
RMS Attenuation Factor = zzzzz	

1. Record the "Magnitude Squared Attenuation Factor" number and record it in the appropriate Value box in Table D5-2.

TABLE D5-2
ATTENUATION FACTORS

MODE	CALIBRATED HARDWARE	PART NUMBER(S)	VALUE	NOMINAL VALUES
BODY OR HEAD	DUMMY LOAD + CABLES ATTEN FACTOR	46-255595G3 46-317724P14 46-255595G4		1000 TO 1200
HEAD	AMP TO WATTMETER ATTEN FACTOR	46-255595G3		1.0 TO 1.2

D-6 Calculation of RF power

Peak voltage should be used in the calculation in order to get an accurate result. It can be converted to power using the following formula as long as certain factors are known and accounted for. The scope correction factor **MUST** be known. So must the *actual* total loss attributed to anything that connects the measuring device to the source. This often includes the accumulated loss associated with the dummy load and any interconnecting cables. Table D6-1 shows the calculation of power if all the attenuating devices in the measurement circuit exhibited perfect loss; that is, the devices added no more or less loss than what they were designed to provide. Table D6-2 shows the same calculation of power but accounts for the measurement-circuit loss values in deriving the true power. Note that the loss has a significant impact on the calculated power.

WARNING!

THE SCOPE CORRECTION FACTOR MUST BE KNOWN FOR THE FORMULAE SHOWN IN TABLES D6-1 AND D6-2. IF IT IS NOT KNOWN, DO NOT USE THIS METHOD. GROSSLY INACCURATE MEASUREMENTS AND POSSIBLE SYSTEM DAMAGE WILL RESULT.

$$\frac{\left(\frac{V_{\text{peak}}}{\text{scope correction factor}} \right)^2}{2 \times Z} \times \text{dummy load and cables attenuation} = \text{RF Power}$$

where "X" and "X" in the above formula signifies multiplication

Assume $Z = 50 \, \Omega$, $V_{\text{peak}} = 27.53$, scope correction factor = 1.00 (scope loss is negligible at 1.0T frequency), dummy load and cables atten. = 1000 (dummy load and cables are all ideal)

$$\frac{\left(\frac{27.53 V_p}{1.00} \right)^2}{100} \times 1000 = 7580 \text{ Watts} = 7.58 \text{ kW}$$

This result assumes a ***theoretically perfect*** situation in which there is no loss. These situations, in common practice, rarely exist!

TABLE D6-1
RF POWER CALCULATION WITH NO LOSS

Now, consider the "real life" type situation in Table D6-2 in which the loss is considered:

$$\frac{\left(\frac{V_{peak}}{\text{scope correction factor}} \right)^2}{2 \times Z} \times \text{dummy load and cables attenuation} = \text{RF Power}$$

where "X" and "X" in the above formula signifies multiplication

Assume $Z = 50 \, \Omega$, $V_{peak} = 27.1$, scope correction factor = 1.00 (scope loss is negligible at 1.0T frequency),
dummy load and cables atten. = 1028

$$\frac{\left(\frac{27.1V_p}{1.00} \right)^2}{100} \times 1028 = 7550 \text{ Watts (very close to 7.58kW.)}$$

Accounting for the loss resulted in an accurate answer. 7550 Watts is as close as we can hope to get to 7580 Watts without using the RF Power Measurement Kit. Note that if none of the loss had been accounted for the error could have been **206 Watts or 2.7%**. As a result, the observer would attempt to adjust the RF power above the 7580 Watt limit.

TABLE D6-2
RF POWER CALCULATION WITH LOSS CONSIDERED

D-7 System Restoration

1. Reconnect original cables to the System Cabinet I/F J1 and J2.
2. Remove the female-BNC to female-BNC (bullet) adapter joining A24 A1 A1 J12 and A24 A1 A1 J2 and reconnect these to the body TNS. J2 (signal input) will connect to the top of the TNS and J12 (signal output) will connect to the bottom.
3. Perform one satisfactory head or body scan.

D-8 Wattmeter Trivia

Why are measurements between a oscilloscope and a Bird Wattmeter found to be different? Repeated GEMS Education Center lab groups that performed the Dummy Load and Cables attenuation process correctly, and obtained an accurate attenuation value, were able to accurately measure RF power as well as with a 400 MHZ scope! Remember, for wattmeter measurements however, the meter reading must be multiplied by the atten factor. Leaving this little factor out is the reason why the above 'note' was made. In body mode the wattmeter is placed after the dummy load. Its atten value must be used to calculate the correct RF value. If it is NOT used, then it should be expected there will be an error of 10% or more!

In head mode, the wattmeter is placed before the dummy load. In this case only the 'cable' factor is used to correct the wattmeter reading. The 'cable' factor is the other atten test you did for the 10 foot cable only. Typically that value is between 1.0 to 1.2. Doesn't sound like much???? Take your calculator out for a moment. Take 2,000 watts (displayed reading) and multiply by an atten factor of 1.0 ; the results are obvious, 2,000 watts. Now take the same displayed 2,000 and multiply it by an atten factor of 1.2 = 2,400 watts! Hopefully this example will show the importance of the attenuation values. You just calculated a 20% error if you were wrong because you didn't do atten test on the 10 foot cable!

D-9 Scope Trivia

A basic 'rule of thumb' in the scope industry is that 'to accurately measure a high frequency signal, the bandwidth of the scope should be at least 5X the frequency being measured'. That means that for a 64 MHZ signal the scope bandwidth would have to be at least 320 MHZ. Along the same line, a 42 MHZ signal would require a bandwidth of at least 210 MHZ. This is why the 400 MHZ scope was chosen several years ago for MR use. Knowing this should alert you as to what might have to be done in order to get correct scope readings from whatever scope you carry today. Who knows what bandwidth scope you'll be working with tomorrow? Will it be of sufficient bandwidth to give accurate measurements? What can you do if it is not covered by the rule of thumb? Read on!

Experimentation with multiple scopes have shown that an error is typically somewhere between 0 and 16%. That's right, some 100 MHZ scopes have no loss! Others may have as much as 15 or 16%. Additionally, it has been found that it is not necessarily the whole scope that has this loss, but the individual channels typically exhibit differences. In other words, channel one might have a 3% loss when channel two could lose 15%. There is no fixed number for all 100 MHZ scopes!

How do you get around these losses? Well, there are two ways. First you could compare a RF waveform on your scope, side by side, with a 400 MHZ scope. Use the same cables etc. and just move the cables from scope to scope. You can then calculate the loss for each channel of the 100 MHZ scope vs. the 400 MHZ scope. Keep the scope intensity low for best accuracy. If and when possible, determine the losses for a 1.5T and 1.0T system. That's for 64 and 42 MHZ. Point of note: there are various losses at 42 MHZ as well as with 64 MHZ. And, typically they differ! Your scope may have a 0 to 15% rolloff at 64 MHZ and 0 to 12% at 42 MHZ. All scopes differ!

Secondly, the next time you have your scope calibrated, just attach a tag asking for the % loss of signal at 42 and 64 MHZ for each channel. You should also state that you want this to be done at the 5 and 2 volts/div setting. This will cover both body and head measurements. At two different re-calibration vendors we use, there are no extra charges incurred for this request.

As mentioned above, the 100 MHZ scope has been shown to be as accurate as a 400 MHZ scopes if the scope rolloff factor is known. Additionally, each channel has it's own unique loss factor.

The biggest error found in GEMS Education Center labs when measuring RF power has proven to be the FE using "typical" dummy load and cables attenuation factors. Basically there are no "typical" factors! All loads and cables are somewhat different. Several may appear to be close in value, but in most cases they are different. If you use a "typical" atten factor in your calculations, then "typically" you will be in error.

Through experimentation at the GEMS Education Center, it has been determined that if you know the attenuation factor for a particular dummy load and cables, those values hold over several years. They do NOT tend to change. This means that if you could do atten test on a dummy load and cables, and keep them as a kit, you don't have to recalibrate them every time you use them. Wow! And possibly a little cheaper than the Power Measurement kit.

If you must, or intend to use a scope to measure body RF power, use atten test on the dummy load and cables.

Knowing the correct attenuation values and the correct scope rolloff values will provide you with measurements as accurate as a 400 MHZ scope and / or the GE Power Measurements kit.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
1	May 19, 1999	Resa Lambert	Updated initial 1.0T procedure. Corrections made as well including new Illustrations for equipment setup. Released on 2160623-1 Rev7 CD-ROM.
2	Sept 2, 1999	Resa Lambert	Updated initial 1.0T procedure to include 1.0T and 1.5T versions, therefore, combining 2 documents into one new document beginning at Rev 0. Corrections made as well including new Illustrations for equipment setup. Added Appendix B, C and D. Added SCOPECAL.xls Tool. Added measurement example/tutorial. Moved Head checks first in the process. Released on 2160623-1 Rev8 CD-ROM.
3	Dec 2, 1999	Resa Lambert	Updated P.M. Tables for SSM vs. RFSC.
4	October 18, 2000	Don Thome'	Created 1.0T, class C specific document. Rewrote RF measurement section and added directions for using RF Calculator Tool. Added process flowchart in Overview. Added Alternate Non-proprietary protocols (Appendix B), Dummy Load and Cable Calibration document (Appendix D).
5	December 20, 2000	Don Thome'	Corrected switch setting in Section 7. Removed reference to the SRF cabinet.
6	April 20, 2001	Don Thome'	Added reference to MONS1.EXE.
7	May 30, 2001	Don Thome'	Simplified power calculation sections for both head and body. Corrected errors throughout document.
8	July 6, 2001	Don Thome'	Added reference to Power Calculator Tool throughout document. Corrected various errors.
9	Oct. 3, 2001	Don Thome'	Corrected error in Head Pulse Width High minimum parameter in Table 5-2.
10	January 4, 2002	Don Thome'	Changed power measurement process so that Vp is measured instead of Vpp to obtain more accurate results. Added information concerning use of MONS1.EXE software.
11	March 1, 2002	Don Thome'	Corrected error in statement after Illustration 4-4 on page 10.