

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- REMOVAL OF FAILED HEAD T/R SWITCH AND PREAMP	2
2- INSTALLATION OF NEW HEAD T/R SWITCH AND PREAMP.....	4
3- FUNCTIONAL CHECKS REQUIRED	5
REVISION HISTORY	6

This procedure applies to 1.5T and 1.0T systems.

1- REMOVAL OF FAILED HEAD T/R SWITCH AND PREAMP

CAUTION

Possible equipment damage. Verify that power is removed from RF cabinet before removing or disconnecting exam room RF components. Verify that the breaker on rear of cabinet is off before starting work. The RF amplifier or TR driver could be damaged by operating into improperly terminated loads.

1. Turn off RF Cabinet power circuit breaker, CB1, on rear of RF Cabinet. Lock out and/or tag out the RF Cabinet.
2. Verify that power to the RF/Pen II Cabinet is off by checking the fans and LEDs on both the RF amplifier and the System Support Module (SSM).
3. Dock and raise table, select new patient, and move the head carriage assembly to the rear of the magnetic bridge.
4. Remove four screws holding carriage cover to carriage assembly. See Illustration L2923A .

STEP 3

REMOVE AND RETAIN FOUR SCREWS AND WASHERS FROM CARRIAGE COVER

STEP 4

SLIDE CARRIAGE COVER FORWARD OVER ACTUATORS, LIFT, AND MOVE BACK TO EXPOSE CARRIAGE HOOK.

STEP 5

REMOVE J6 AND THE PREAMP OUTPUT CABLE.

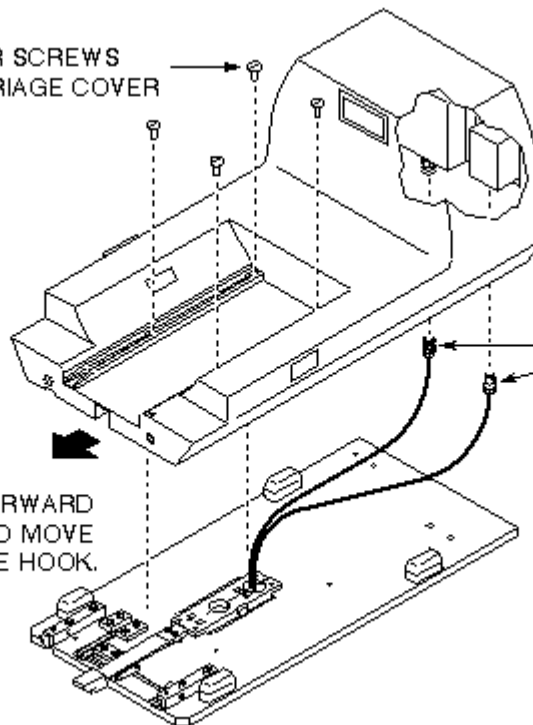


ILLUSTRATION L2923A
REMOVAL OF CARRIAGE COVER

L2923A

5. Remove cover from assembly.

6. Remove cable from J6 of the T/R Switch and the cable from the preamp output.
7. Turn carriage on its back. See Illustration L2924A .

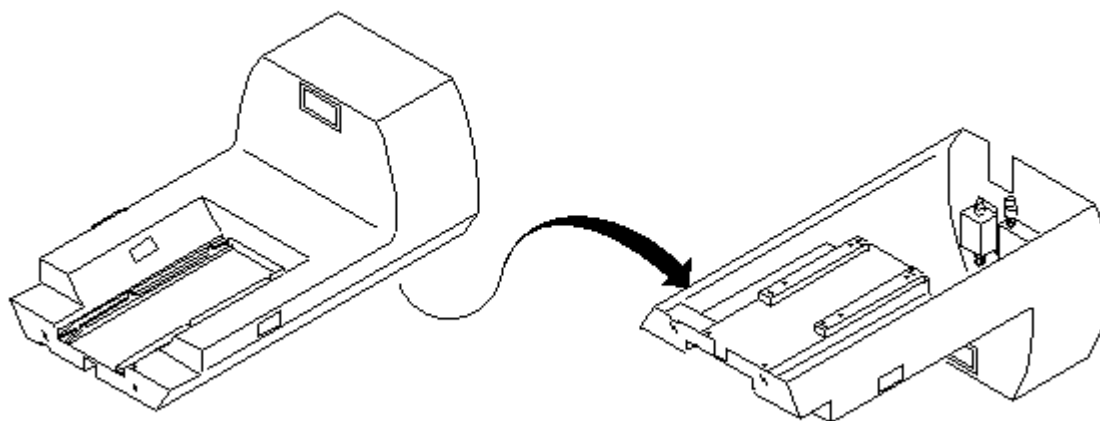


ILLUSTRATION L2924A
INSIDE CARRIAGE COVER

8. Disconnect cable between T/R Switch and preamp input. See Illustration L2925A .

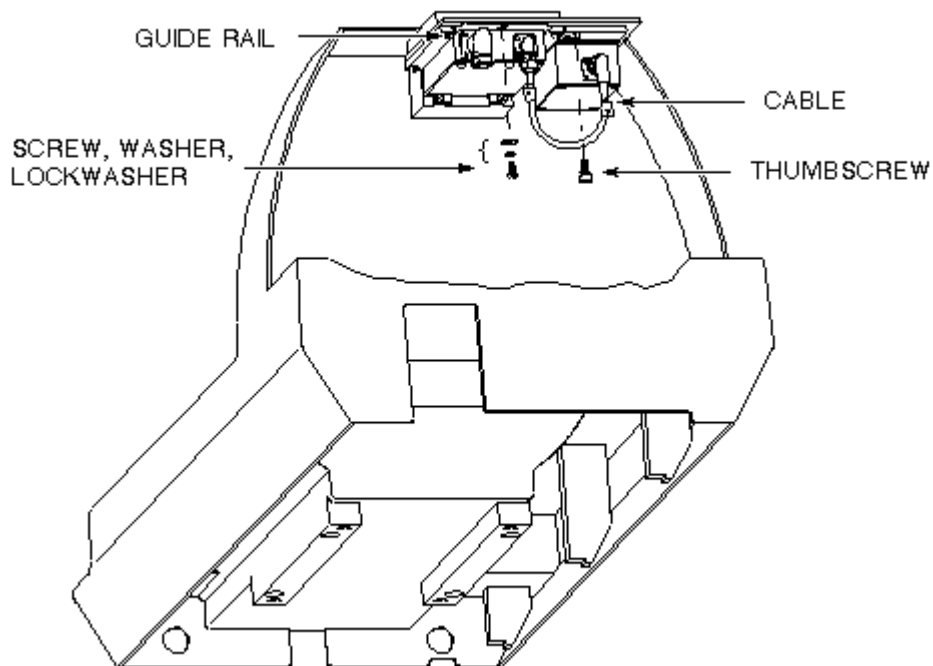


ILLUSTRATION L2925A
HEAD T/R SWITCH AND PREAMP

9. Remove the T/R Switch (screw, washer and lockwasher) and/or preamp (thumb screw).

2- INSTALLATION OF NEW HEAD T/R SWITCH AND PREAMP

CAUTION

Possible equipment damage. Verify that power is removed from RF cabinet before removing or disconnecting exam room RF components. Verify that the breaker on rear of cabinet is off before starting work. The RF amplifier or TR driver could be damaged by operating into improperly terminated loads.

1. Install the T/R Switch (screw, washer and lockwasher) and/or preamp (thumb screw). See Illustration L2925A .

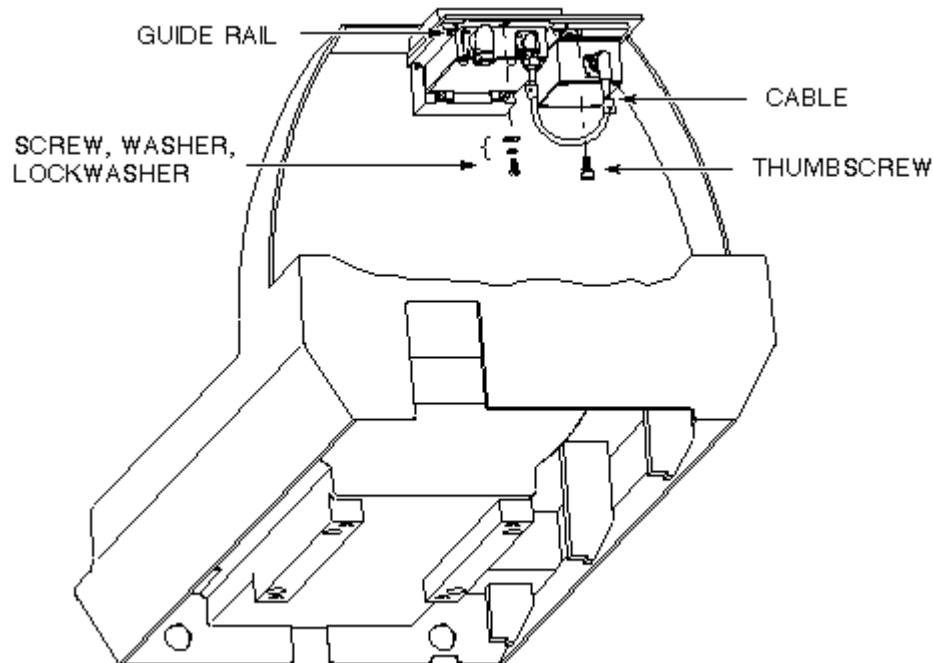


ILLUSTRATION L2925A
HEAD T/R SWITCH AND PREAMP

2. With T/R Switch assembly fully inserted, tighten screws on the guide rails. Tighten preamp thumbscrew.
3. Connect cable between T/R Switch and preamp input.

See Illustration L2927A for steps 4-6.

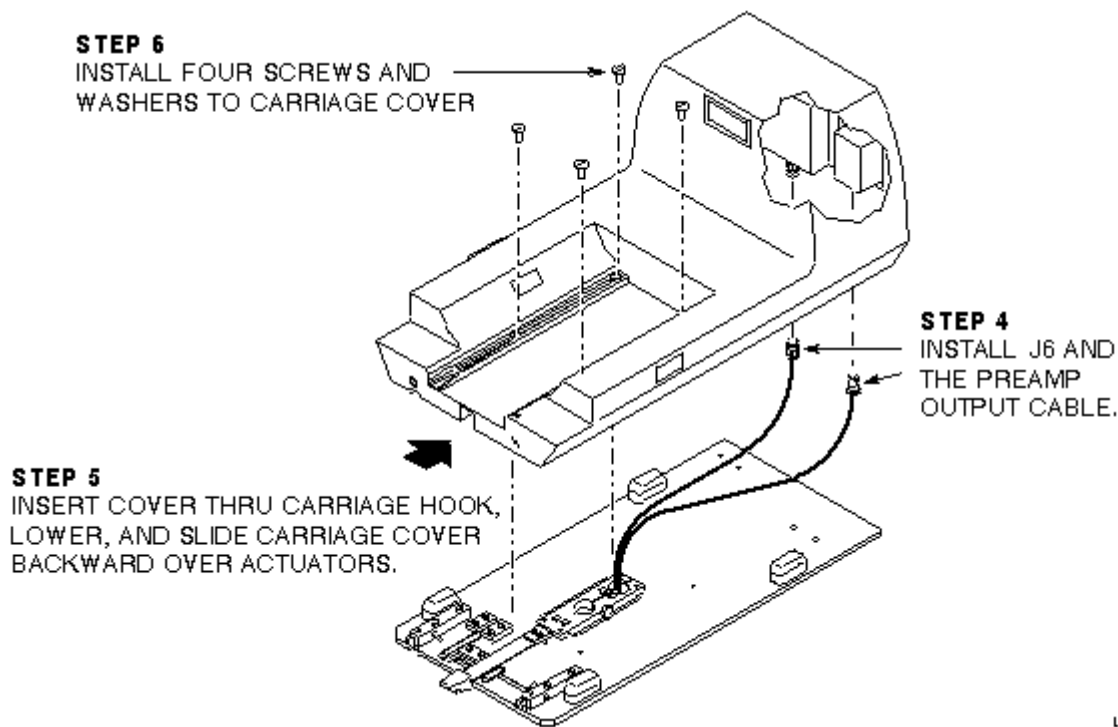


ILLUSTRATION L2927A
INSTALLATION OF CARRIAGE COVER

4. Connect cable to J6 of the T/R Switch and the cable to the preamp output.
5. Install cover to assembly.
6. Install 4 screws to hold carriage cover to carriage assembly.
7. Remove lock out and/or tag out from the RF Cabinet. Turn on RF Cabinet power circuit breaker, CB1, on rear of cabinet.

3- FUNCTIONAL CHECKS REQUIRED

1. Refer to Functional Checks tab, Procedure DD/TR Driver And Pin Diodes.
2. Check Head SNR. Refer to procedure for Signal To Noise Check.

(An alternative proprietary procedure is available for GE use or to customers with an Advanced Service Package Limited License. Refer to procedure for TLT and perform a Head SNR scan.)

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	Aug 19, 1998	R. Hawthorne	Initial conversion to Word