

TABLE OF CONTENTS

TABLE OF CONTENTS	1
1- DICOM PRINT CONFIGURATION	2
1-1 Description	2
1-2 Set-up Procedure	2
1-3 Advanced DICOM Camera Attributes (Release 8.3 or later only)	7
2- DICOM NETWORK CONFIGURATION	9
2-1 Description	9
2-1 Set-up Procedure	9
3- DICOM ARCHIVING CONFIGURATION	12
3-1 Description	12
3-2 Set-up Procedure	12
REVISION HISTORY	14

1- DICOM PRINT CONFIGURATION

Note:

When upgrading from Release 8.2.5 to 8.3, some vendor dicom print equipment will need to be reconfigured to accept the new DICOM print AETitle. The old AETitle was ***Print_SCU***; The new AETitle is ***<hostname>_DCP*** where *<hostname>* is the hostname of the computer (e.g. a hostname of *signa* has an AETitle of *signa_DCP*). This information must be given to the camera vendor to input the new AETitle name for cameras like the Kodak/Imation 9410 that require the AETitle be pre-configured.


1-1 Description

The DICOM Print feature allows configuration of the filming function to accommodate filming to any paper or film printer that is networked to the workstation through a DICOM printer.

1-2 Set-up Procedure


Access to DICOM print set-up is done via the Display desktop. To configure the DICOM printing:

1. Enter the Display Desktop. The Browser window will open.
2. Select **[Film Composer]** on the Browser window. The Film Composer pop-up window will appear (refer to Illustration 1-1).

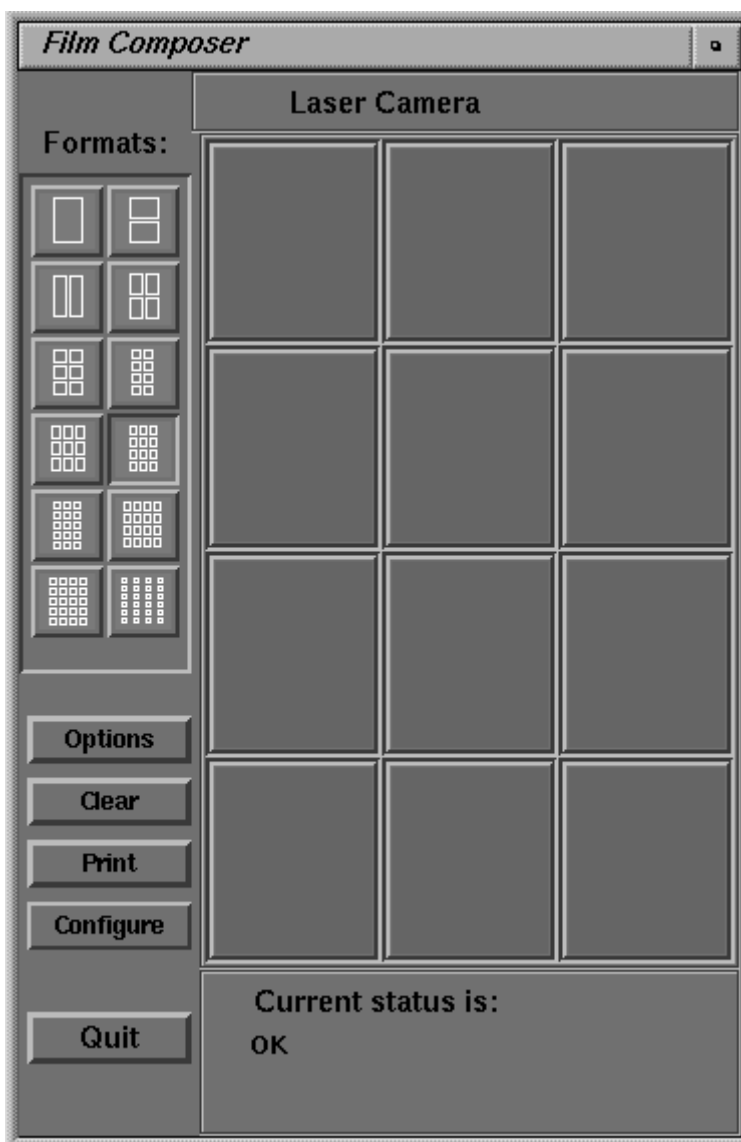

FILM COMPOSER WINDOW
ILLUSTRATION 1-1

3. Select **[Configure]**, the Remote printer configuration pop-up window appears (refer to Illustration 1-2). Table 1-1 lists the descriptions for each of the buttons on the pop-up window.

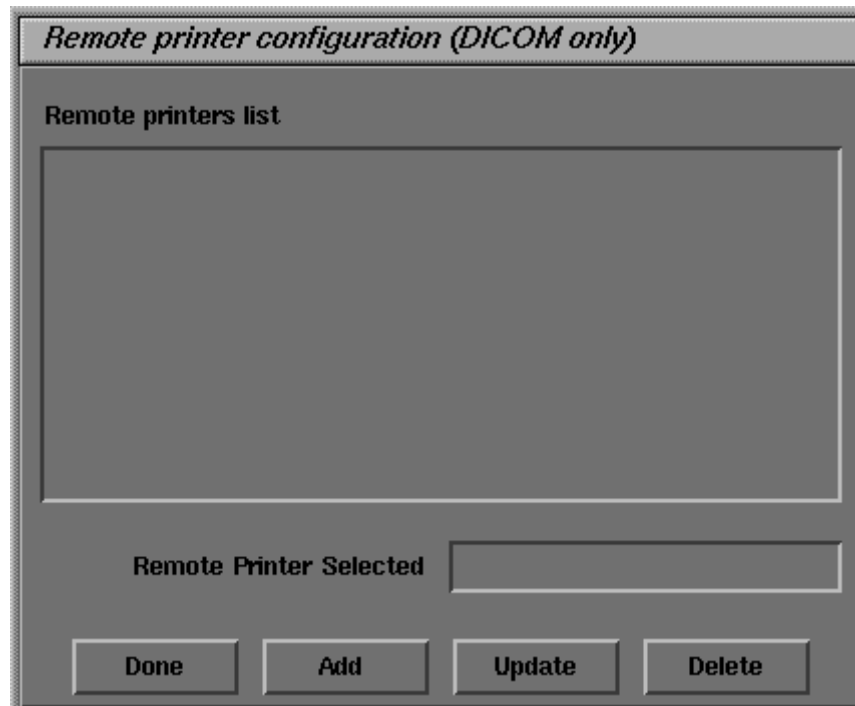

REMOTE PRINTER CONFIGURATION
ILLUSTRATION 1-2

TABLE 1-1
REMOTE PRINTER CONFIGURATION BUTTONS

Button Name	Description
[Add]	Brings up the Camera Type popup, select type to bring up Configuration Parameters pop-up window to allow printer configuration.
[Delete]	Erases the highlighted printer from the list
[Update]	Opens the highlighted printer configuration window and allows modifications to the existing parameters.
[Done]	Accepts any changes and additions made and closes the window.


4. Select **[Add]**, then **[Dicom]** Camera type. See Illustration 1-3. This brings up the Camera Configuration Dialog Box with predefined camera types. See Illustration 1-4.



SELECT CAMERA TYPE
ILLUSTRATION 1-3

Note

The DASM button on the “Type of Laser Camera” window is only activated if the system has a DASM module and the FE has installed the DASM via the Guided Install Tool.

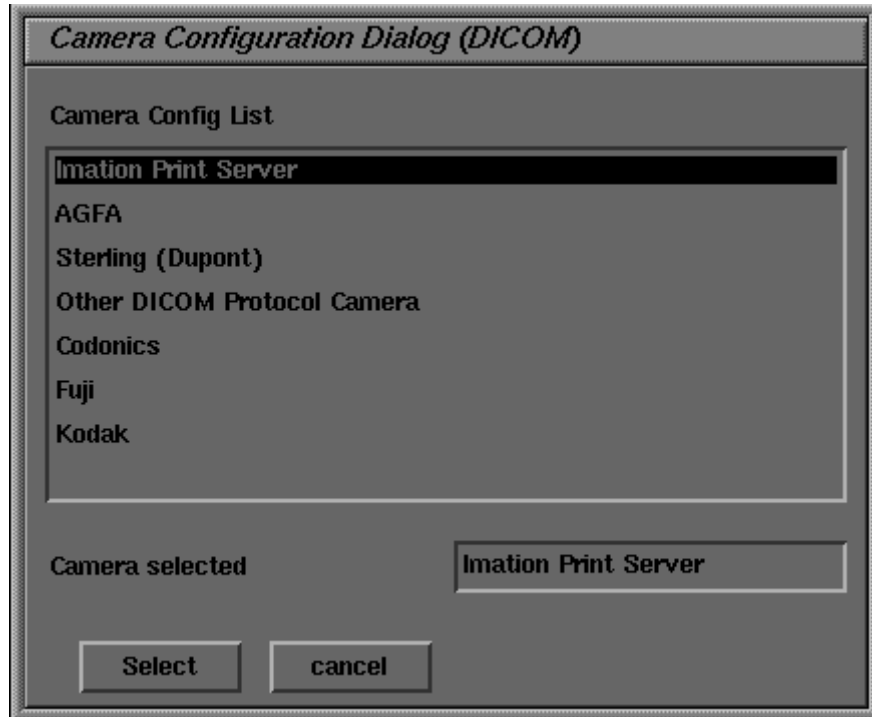

CAMERA CONFIGURATION DIALOG BOX
ILLUSTRATION 1-4

4. Select the camera configuration to match the camera type. The DICOM printer configuration pop-up window shown in Illustration 1-5 appears. This window is used to configure a printer and set filming parameters.

DICOM printer configuration parameters

Device Name:

Host Name:

IP address:

Application title:

TCP listen port:

Comments:

Medium Type:

Destination:

Film orientation:

Magnification type:

Standard film formats:

F_1x1	F_2x1	F_2x2	F_3x2
On	On	On	On
F_3x3	F_4x2	F_4x3	F_4x4
On	Off	On	On
F_5x3	F_5x4	F_6x4	
On	On	On	

THESE PARAMETERS ADDED AT RELEASE 8.3

DICOM PRINT CONFIGURATION
ILLUSTRATION 1-5



The default parameter settings have been selected to match the capability of the different cameras supported by the software. Keep in mind that the interface is not programmed to keep you from making a bad or incorrect selection. For example, an AGFA Camera does not support a 4 x 2 format, hence this format is defaulted to OFF. However, you can turn this illegal format ON and the GUI will show it as on, but the camera will still not be able to print film in that format.

- **Device Name** – (For Release 8.3 or greater). Name for the printer (name preferred by the customer)
- **Host Name** – Network Host name for the printer (supplied by the customer)
- **IP Address** – (For Release 8.3 or greater). Network IP Address that corresponds to the Host Name for printer (supplied by the customer)
- **Application title** - The DICOM name of the printer (provided on the DICOM conformance statement for the camera; supplied by the camera FE)
- **TCP listen port** - Port which the DICOM printer accepts print requests (provided on the DICOM conformance statement for the camera; supplied by camera FE)
- **Comments** - customer reference only
- **Medium Type** - {Paper, Blue Film, Clear Film}. Type of media.
- **Destination** - {Magazine, Processor}. Destination of the printed media.
- **Film orientation** – {Portrait, Landscape (currently not selectable)}. Orientation of the film/paper.
- **Magnification type** – {Replicate, Bilinear, None, (or Cubic for Release 8.3 or greater)}. This is a customer preference or as specified by the camera FE.
- **Standard film formats** - Select [On] for each format option the customer wants listed on the Film Composer when the printer is selected. Note: the default settings for film formats are “On” if supported by that camera.
- **Advanced** – (For Release 8.3 or later) Brings up Camera Attribute Dialog Box. See Section 1-3 for description of Advanced settings.

4. Select **[OK]** to accept the selections and close the window.
5. Beginning with Release 8.3, SaveINFO will save the DICOM print configuration settings. To perform a SaveINFO:
 - a) Select the Service Desktop, select **[Install]**.
 - b) Enter **operator** for password
 - c) After the Install GUI is displayed, select **Save/Restore** tab.
 - d) Insert InfoMOD and select **SaveInfo**.

1-3 Advanced DICOM Camera Attributes (Release 8.3 or later only)

Beginning with Release 8.3, the ability to customize the default printer settings has been added to the software. FE's working closely with a camera manufacturer representative and the customer can modify the default settings to adjust film output to individual customer preferences. See Illustration 1-6.



The default attribute settings have been selected to match the capability of the different Cameras supported by the software. Keep in mind that the interface is not programmed to keep you from making a bad or incorrect selection. For example, a camera that does not support “Smoothing Type” will have this setting defaulted OFF. However, you can turn this illegal type ON and the GUI will show it as on, but the camera will still not be able to apply Smoothing to the images.

A screenshot of the "Camera Attribute Dialog" window. The title bar says "Camera Attribute Dialog". Inside, under "Image Quality Parameters", there are several settings: "Configuration" (empty text box), "Maximum Density" (text box with "320"), "Minimum Density" (text box with "32"), "Smoothing Type" (radio buttons for "On" and "Off", with "Off" selected), "Empty Density" (radio buttons for "Black" and "White", with "Black" selected), "Border Density" (radio buttons for "Black" and "White", with "Black" selected), "Film Size" (radio buttons for "<empty>" and "No", with "<empty>" selected), "Trim" (radio buttons for "No" and "Yes", with "No" selected), and "Priority" (radio buttons for "High" and "Low", with "High" selected). At the bottom are "Done" and "cancel" buttons.

CAMERA ATTRIBUTE DIALOG BOX (RELEASE 8.3 ONLY)
ILLUSTRATION 1-6

- **Configuration** – *default is to leave this empty.* This field provides a means to enter very specific camera configurations strings provided by the camera manufacturer to modify the specific camera settings.
- **Maximum Density** – *default 320.* These values are site specific and/or camera installation dependent. Some information may be specified in the camera's conformance statement. More often, the values are defined by the camera installer and require knowledge/information stored on the camera. The site may have specific requirements or settings.
- **Minimum Density** – *default On = 32.* {On, Off}. See usage as described above for Maximum Density.

- **Smoothing Type** - *default is to leave this empty.* Requires Magnification Type = Cubic on Dicom Print Configuration Parameters screen. The site may have specific requirements or settings.
- **Empty Density** – *default is Black or Off if not supported by camera.* {Off, Black, White}. “Off” means this element and value is not sent, as some cameras may not support this parameter and sending the element causes a camera error. Consult the camera's conformance statement for support of this parameter.
- **Border Density** – *default is Black or Off if not supported by camera.* {Off, Black, White}. “Off” means this element and value is not sent, as some cameras may not support this parameter and sending the element causes a camera error. Consult the camera's conformance statement for support of this parameter.
- **Film Size** – *default is empty {empty, 14INX17IN}.* Recommended selection is <empty>. Some sites may require specifying the film size as the camera has multiple film supplies with different film sizes.
- **Trim** – *default is No or Off if not supported by camera.* {Off, Yes, No}. “Off” means this element and value is not sent, as some cameras may not support this parameter and sending the element causes a camera error. Consult the camera's conformance statement for support of this parameter.
- **Priority** – *default is High {High, Medium, Low}.* This is the relative priority given to the task of printing vs. all other processes running on the scanner.

Beginning with Release 8.3, SaveINFO will save the DICOM print configuration settings. To perform a SaveINFO, refer to step 5 in previous section.

2- DICOM NETWORK CONFIGURATION


2-1 Description

The DICOM network allows automatic download of exams or series to a selected node(s).

2-1 Set-up Procedure

Access to DICOM networking set-up is done via the Display desktop. To configure the DICOM network:

1. Enter the Display Desktop. The Browser window will open.
2. From the tool bar, select the “**Network**” pull-down menu using the left mouse button. Then select “**Remote Host**”. The following window will appear (see Illustration 2-1).

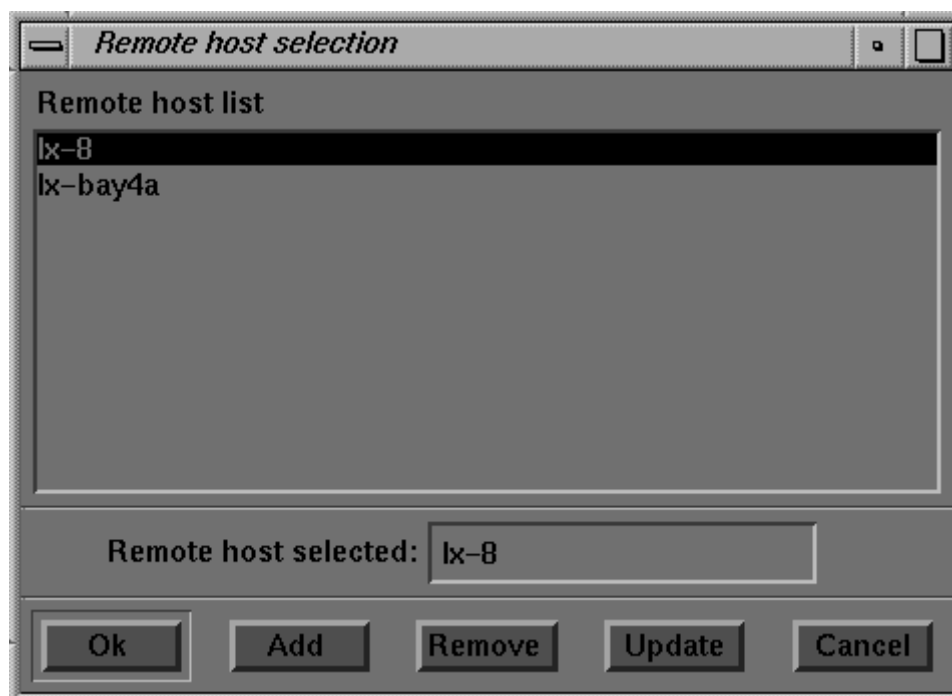

REMOTE HOST SELECTION
ILLUSTRATION 2-1

Refer to Table 2-1 for a list of the button names and their functions.

TABLE 2-1
REMOTE HOST CONFIGURATION BUTTONS

Button Name	Description
[OK]	Accepts changes and closes the window
[Add]	Opens the DICOM Printer Configuration Parameters pop-up window to allow printer configuration.
[Delete]	Erases the highlighted printer from the list
[Update]	Opens the highlighted DICOM printer configuration window and allows modifications to the existing parameters.
[Cancel]	Cancels the window without making changes

- To add a remote host, select **[Add]**. Illustration 2-2 shows the pop-up window that will appear

Remote host parameters

Host name

Network address

Network protocol

Advantage NET
DICOM

Port number

4005

AE Title

Comment

Archive Node: ☒ Auto ☐ Yes ☐ No

Save Cancel

REMOTE HOST PARAMETERS
ILLUSTRATION 2-2

- Host Name - Enter the host name (supplied by the customer)
- Network Address - Address of Host Name (supplied by customer)
- Network protocol - Select DICOM or Advantage Net (if this will access other GE equipment, select Advantage Net)
- Port number - supplied by customer
- AE Title – Application Entity for DICOM standard. It provided by the remote device.
- Comment - for user reference
- Archive Node - If the user would like a remote host to be used to remotely archive, then select [Yes], the default is set as [No].

4. Select **[Save]** to accept the changes and close the window. Then press **[OK]** in the previous window.

3- DICOM ARCHIVING CONFIGURATION

3-1 Description

The hardware will automatically set-up what devices are connected to the system for archiving purposes. When the user selects archive, the selected images/exams/series will be automatically archived to the default selection.


3-2 Set-up Procedure

The archive device selection is accessed via the Display desktop. To select the default media for archiving:

1. Enter the Display Desktop. The Browser window will open.
2. From the tool bar, select the “**Archive**” pull-down menu using the left mouse button. Then select “**Remote Host**”. The following window will appear (see Illustration 3-1).
3. Select the device that the customer wants selected as the default and press the [**Set Default Device**], then press [**OK**].

Note:

The example shown in Illustration 3-1 shows only one device found by the hardware. The site may have more than one.

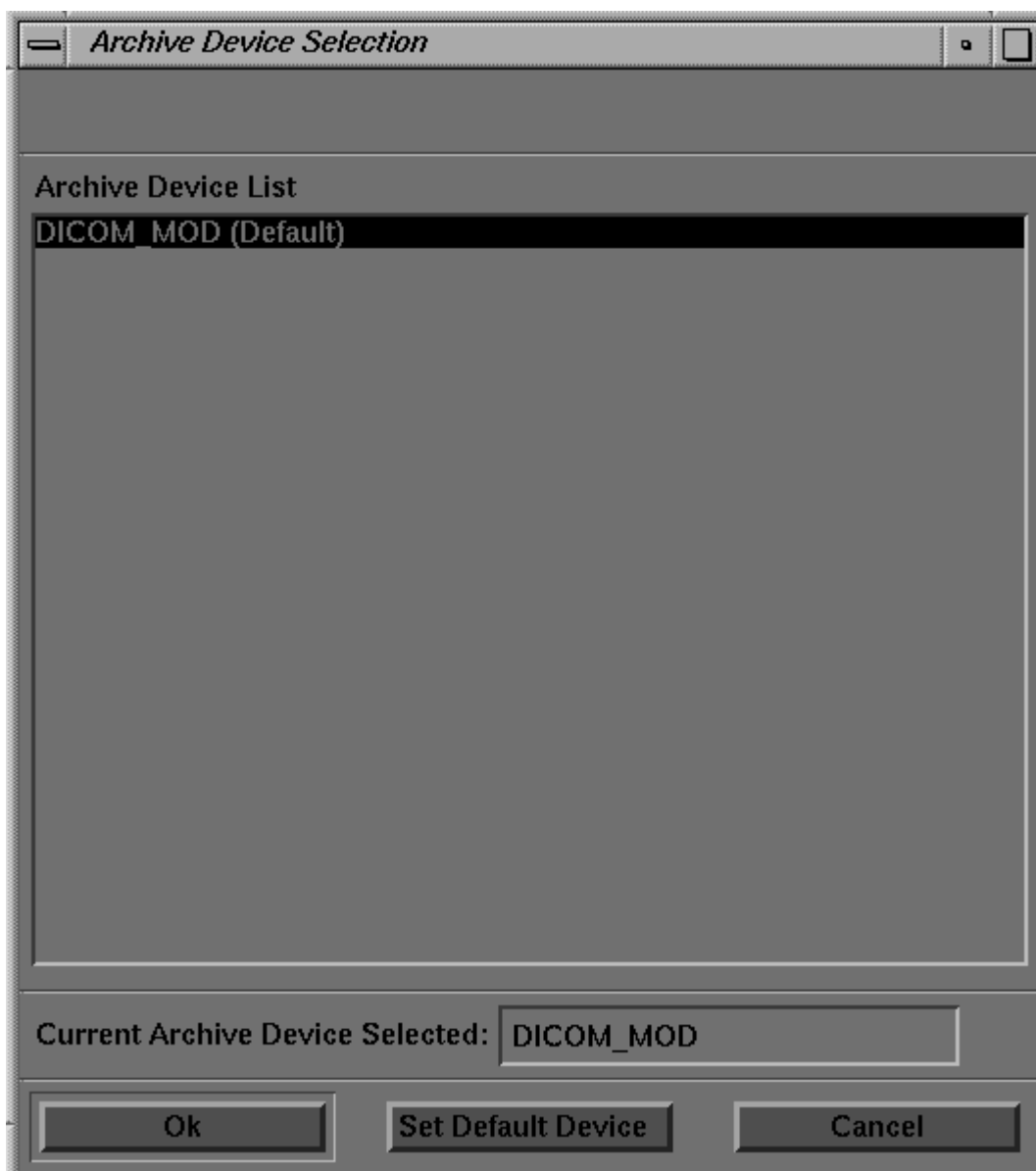

ARCHIVE DEVICE SELECTION
ILLUSTRATION 3-1

3. Select the device that the customer wants selected as the default and press the **[Set Default Device]**, then press **[OK]**.

REVISION HISTORY

REV	DATE	AUTHOR	PRIMARY REASONS FOR CHANGE
0	May 5, 1998	K. L-P	Initial release
1	Mar. 2, 2000	J. Wolak	Updated for Release 8.3 new features and impact