

TABLE OF CONTENTS

table of contents ..... 1

1- Introduction ..... 2

    1-1 Procedure..... 2

revision history ..... 7

## 1- INTRODUCTION

Image Snap is a tool used to aid in troubleshooting image quality issues. An image can be captured and stored in a data file for future transfer and review to the Automated Support Center (ASC). The following will be captured into a file:


- Image data files
- Raw data files (if they exist)
- System operating information
- Comments entered by operator

### Note

Insite is required for auto transfer to occur. However, the files will be saved if Insite is not available and/or there is a failure during auto transfer.

### 1-1 Procedure

Image Snap is accessed by selecting **[ImageSnap]** from the Service Desktop. (see Illustration 1-1)

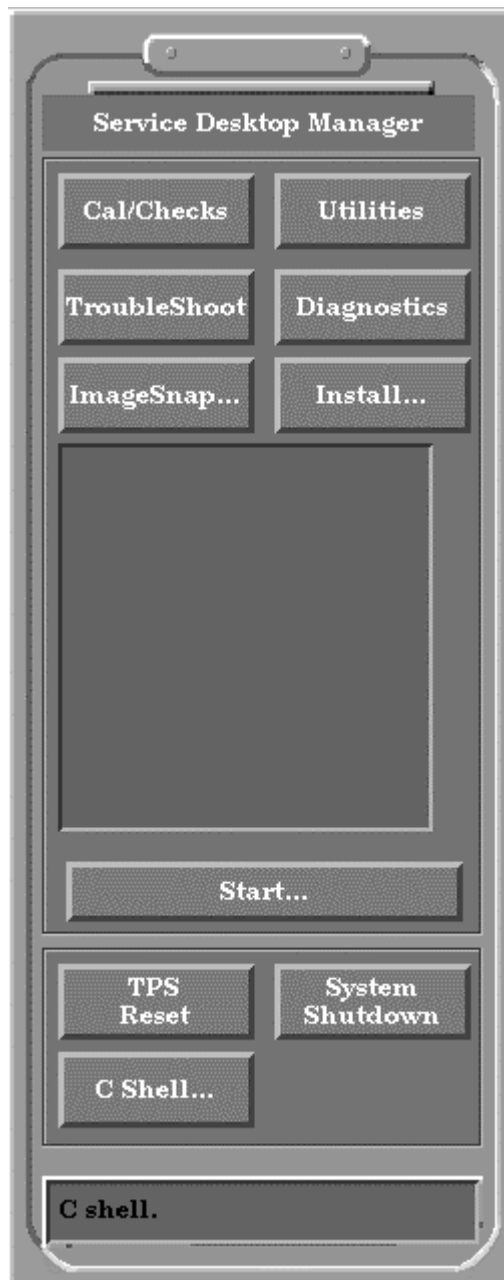

**SELECTING IMAGE SNAP**  
ILLUSTRATION 1-1

After selecting **[ImageSnap]**, a screen opens with available exams and series on the system disk (refer to Illustration 1-2). The most current exam and series is automatically selected. If another exam or series is desired, click on that exam or series. If Raw Data exists, the Yes box will be highlighted. Place cursor in the Comment field to enter comments; files will not be saved without an entry in the Comment field. Only one image file can be selected. A dialog box will be displayed if more than 1 image file is selected.

**Note**

The **[Save]** button will not be selectable until an image is highlighted.

The screenshot shows a window titled "Image Snap Shot". It contains three columns: "Exam:", "Series:", and "Image:". The "Exam:" column has a list of exams with dates and names, where "99 07/09/98 GP FLEX" is selected. The "Series:" column has "1" selected. The "Image:" column has "1" selected. Below these columns is a section for "Raw Data Captured:" with "Yes" and "No" checkboxes. At the bottom is a "Comments from User:" text area and buttons for "Save", "File...", and "Exit".

| Exam:                     | Series: | Image: |
|---------------------------|---------|--------|
| 99 07/09/98 GP FLEX       | 1       | 1      |
| 98 07/09/98 qa quad cspir |         |        |
| 97 07/09/98 CTL REVB      |         |        |
| 96 07/08/98 test          |         |        |
| 50527 07/07/98 CTL SNR    |         |        |
| 50526 07/07/98 CTL SNR    |         |        |
| 50525 07/07/98 CTL SNR    |         |        |
| 50524 07/07/98 CTL SNR    |         |        |
| 50523 07/07/98 CTL SNR    |         |        |
| 50522 07/07/98 CTL SNR    |         |        |
| 50521 07/07/98 CTL SNR    |         |        |
| 95 07/07/98 charlie       |         |        |
| 94 07/07/98 donna         |         |        |
| 93 07/07/98 rhonda        |         |        |
| 92 07/07/98 rhonda        |         |        |
| 91 07/07/98 qa            |         |        |
| 90 07/06/98 andy          |         |        |
| 89 07/06/98 dan           |         |        |
| 88 07/06/98 dan           |         |        |

Raw Data Captured:  
☐ Yes ☐ No

Comments from User:

Save File... Exit

Please select data you want to save

EXAM AND SERIES SELECTION  
ILLUSTRATION 1-2

Verify that the selected file(s) are correct and select **[Save]**. File(s) will be compressed and stored in the directory /usr/g/service/log. A call will be generated to the ASC to notify them (as an email) that an image snap file exists. Select **[Exit]** to return to normal operation. A dialog box will be displayed if more than 10 files exist or if the memory allocation is insufficient to save the selected file.

Select **[File ...]** to verify that the file has been saved (see Illustration 1-3). The files can be deleted from this screen if necessary.



**STORED FILES SCREEN**  
ILLUSTRATION 1-3

If the call to the ASC fails, a message appears on the screen notifying the operator that a call to GE Cares must be done (refer to Illustration 1-4).

The screenshot shows the 'Image Snap Shot' window. It has three columns: 'Exam:', 'Series:', and 'Image:'. The 'Exam:' column contains a list of exams with dates and names. The 'Series:' and 'Image:' columns each contain the number '1'. Below the columns is a 'Raw Data Captured:' section with 'Yes' and 'No' checkboxes. Below that is a 'Comments from User:' text area containing the word 'Test'. At the bottom are 'Save', 'File...', and 'Exit' buttons. A dark status bar at the very bottom contains the text 'Notifying GE failed. Please contact GE CARES' and 'Data is saved at: /usr/g/service/log/bib.lx-bay4a.9807091552'. An arrow points from the text 'Failed message appears here.' to the 'Notifying GE failed...' message.

| Exam:                     | Series: | Image: |
|---------------------------|---------|--------|
| 99 07/09/98 GP FLEX       | 1       | 1      |
| 98 07/09/98 qa quad cspir |         |        |
| 97 07/09/98 CTL REVB      |         |        |
| 96 07/08/98 test          |         |        |
| 50527 07/07/98 CTL SNR    |         |        |
| 50526 07/07/98 CTL SNR    |         |        |
| 50525 07/07/98 CTL SNR    |         |        |
| 50524 07/07/98 CTL SNR    |         |        |
| 50523 07/07/98 CTL SNR    |         |        |
| 50522 07/07/98 CTL SNR    |         |        |
| 50521 07/07/98 CTL SNR    |         |        |
| 95 07/07/98 charlie       |         |        |
| 94 07/07/98 donna         |         |        |
| 93 07/07/98 rhonda        |         |        |
| 92 07/07/98 rhonda        |         |        |
| 91 07/07/98 qa            |         |        |
| 90 07/06/98 andy          |         |        |
| 89 07/06/98 dan           |         |        |
| 88 07/06/98 dan           |         |        |

Raw Data Captured:  
☐ Yes ☐ No

Comments from User:  
Test

Save File... Exit

Notifying GE failed. Please contact GE CARES  
Data is saved at: /usr/g/service/log/bib.lx-bay4a.9807091552

Failed message appears here.

**FAILED MESSAGE EXAMPLE**  
ILLUSTRATION 1-4

## REVISION HISTORY

| REV | DATE      | AUTHOR     | PRIMARY REASONS FOR CHANGE                                                                                            |
|-----|-----------|------------|-----------------------------------------------------------------------------------------------------------------------|
| A   | 12/4/97   | W. Pasciak | Initial Release                                                                                                       |
| 0   | Not known | Not known  | Revision History was not updated for Rev 0 version.                                                                   |
| 1   | 7/9/98    | M. Keber   | Corrected file name in title to indicate it is a class M (not "C") document.<br>Corrected screen illustrations (KL-P) |
| 2   | 6/17/99   | G. Boerner | Minor edits based on 8.3 tools validation                                                                             |
|     |           |            |                                                                                                                       |