
GE Healthcare

Cradle Release Tension and Front Roller Guide Adjustments

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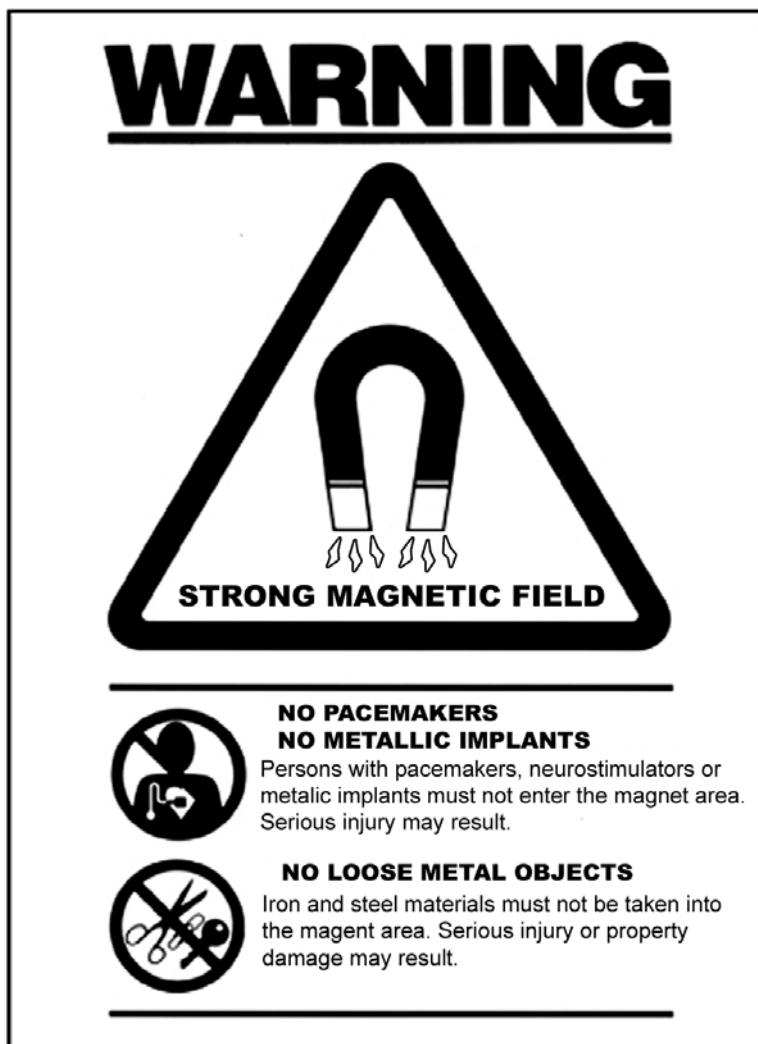WARNING**(EN)**

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR, OR PATIENT FROM ELECTRIC SHOCK, OR FROM MECHANICAL OR OTHER HAZARDS.

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This page left intentionally blank.

Revision History

REV	Date	Author	Primary Reason For Change
1	Oct. 17, 2007	K.Keshena	Initial Release for PSR 13108468, Emergency Cradle Release.

This page left intentionally blank.

Table of Contents

1	Cradle Release Tension Adjustment and Front Guide Roller Height Adjustment	17
1.4.1	Cradle Release From Patient Transport Tension Adjustment	17
1.4.1.1	Cradle Position During Release Mechanism Adjustment	17
1.4.1.2	Initial Checks To Adjust The Cradle Release Mechanism	20
1.4.1.3	Tension Adjustment of Cradle Release Mechanism	22
1.4.2	Front Guide Rollers Height Adjustment	26

This page left intentionally blank.

1 Cradle Release Tension Adjustment and Front Guide Roller Height Adjustment

1.1 Personnel Requirements

Required Persons	Preliminary Reqs	Procedure	Finalization
2			

1.2 Overview

This procedure provides the instructions for the Cradle Release from the Transport Tension adjustment and the front guide roller height adjustment in the patient table.

1.3 Preliminary Requirements

1.3.1 Tools and Test Equipment

Item	Qty	Effectivity	Part#	Manufacturer
Standard Screwdriver	1	-	-	-
3/32 inch Allen Key	1	-	-	-
5/32 inch Allen Key	1	-	-	-
Socket Extension, Nut Driver	1	-	-	-

1.3.2 Consumables

Item	Qty	Effectivity	Part#	Manufacturer
Shims	Amount Needed	-	46-317181P1	-

1.3.3 Safety

WARNING

PERSONAL INJURY.
TABLE HEIGHT IS ADJUSTABLE.
MAKE SURE THE SAFETY LOCK BAR IS IN PLACE BEFORE SERVICING THE PATIENT TABLE. DO NOT LOAD THE TABLE WHILE SERVICING.

1.4 Procedure

1.4.1 Cradle Release From Patient Transport Tension Adjustment

1.4.1.1 Cradle Position During Release Mechanism Adjustment

1. Make sure cradle is positioned on Table.

2. Remove Table from Scan Room to perform rest of service procedure.
3. Press the cradle release pin and front Delrin stopper. See [Illustration 1](#) .

Illustration 1: Location of Cradle Release Pin and Delrin Stopper


4. Pull the cradle out. The cradle will come about 50 mm away from the tabletop front. This is done to make sure the cradle is clear of the flipper actuation. See [Illustration 2](#) .


Illustration 2: Cradle Cleared of Flipper Actuation

5. Press the Up Limit Indicator Bumper to release the central pin. See [Illustration 3](#) .

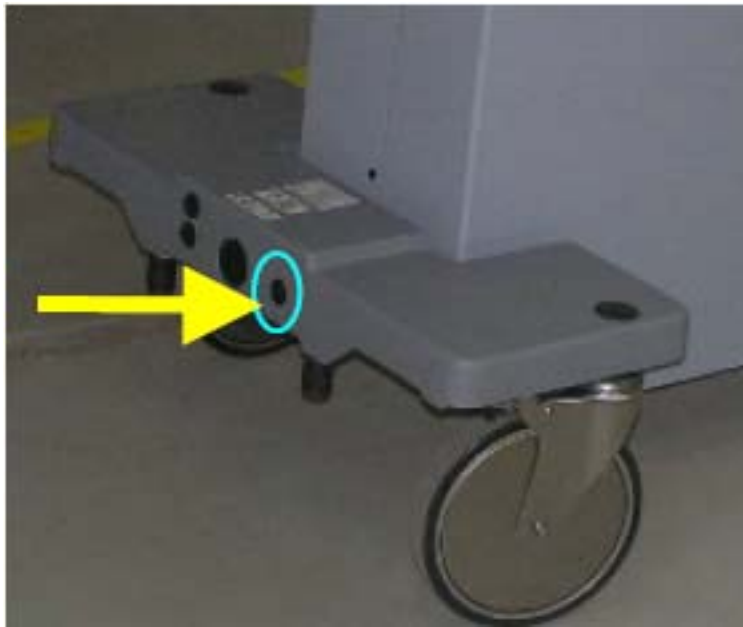Illustration 3: Up Limit Indicator Bumper

6. Pull cradle out and remove from the table. See [Illustration 4](#) .

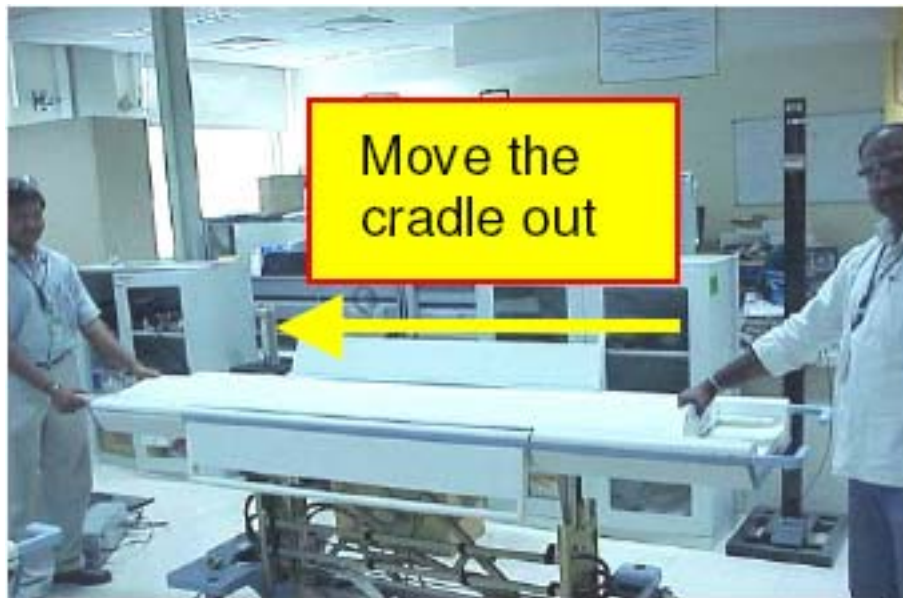Illustration 4: Pull Cradle Off of Table

7. Once cradle has been removed from the table, place the cradle on a surface where further adjustments can be performed. See [Illustration 5](#) .

Illustration 5: Cradle Put Aside

1.4.1.2 Initial Checks To Adjust The Cradle Release Mechanism


NOTICE

Two persons are required to conduct this check

1. Twist Handle fully to extend release bar/ramp as shown in [Illustration 6](#) and hold it at the extended position.

Illustration 6: Emergency Cradle Release

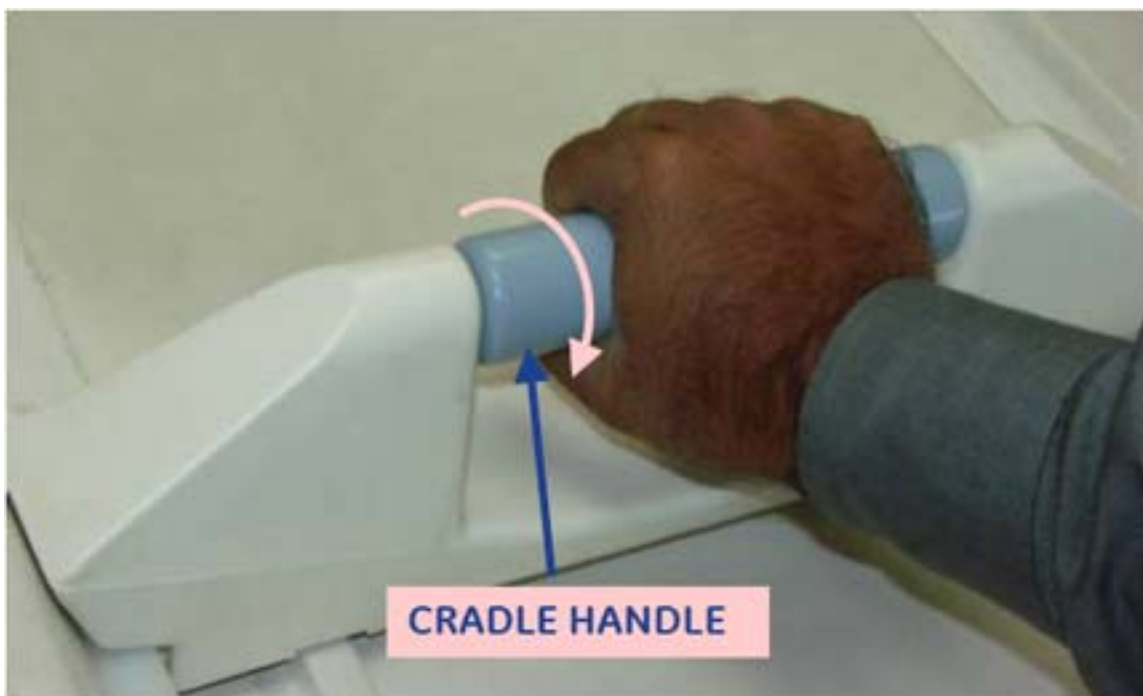2. Push the "Release Bar/Ramp" as shown in [Illustration 8](#) using a socket extension, nut driver, as shown in the illustration, or similar tool with approx. 10-20 lbs force (Force Gauge, 5142800 or 5142800-2). Observe the movement of the release ramp: a properly adjusted release ramp will exhibit some spring and when released will return to its original position.

Illustration 7: Release Bar/Ramp

Illustration 8: Pushing the Release Bar/Ramp

3. If the "Release Bar/Ramp" fails to return to its original position and moves more than Approx. 1.5mm (.06 in.), see [Illustration 9](#) , then adjust the release mechanism tension per instructions in [Section 1.4.1.3](#) .

Illustration 9: Failure Example

1.4.1.3 Tension Adjustment of Cradle Release Mechanism

1. Turn cradle over to access the kevlar cord. See [Illustration 10](#) .


Illustration 10: Location of Kevlar Cord

2. Remove the kevlar cord knot as shown in [Illustration 11](#) . The 3/32 inch Allen Key may be helpful in removing the knot.

Illustration 11: Removing Kevlar Cord Knot

3. Using the 3/32 inch Allen Key, loosen the grub screw as shown in [Illustration 12](#) .

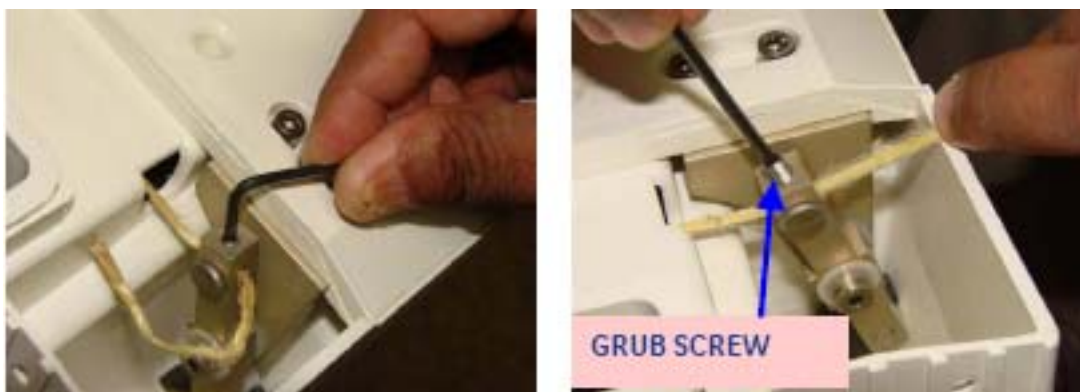Illustration 12: Loosening the Grub Screw

4. Pull the kevlar cord in the direction of the blue arrow shown in [Illustration 13](#) , in order to apply needed tension so the Release Bar/Ramp butts against the cradle front section. See [Illustration 14](#) .

Illustration 13: Securing Tension On The Kevlar Cord

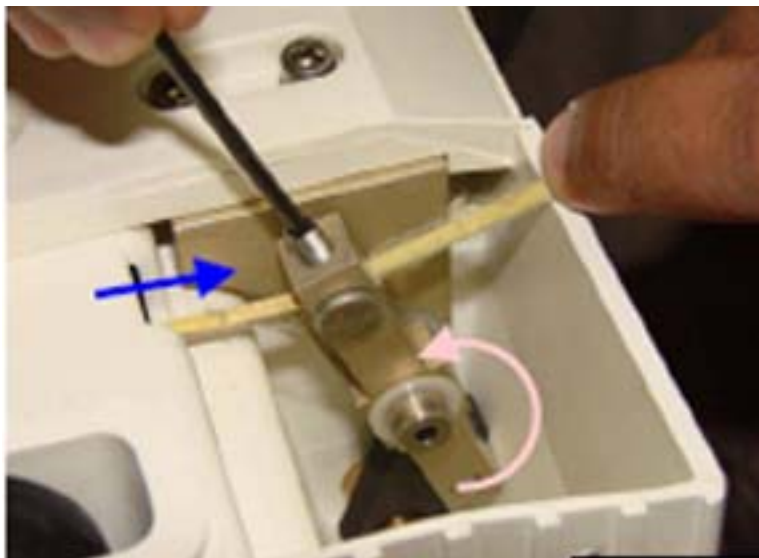
Illustration 14: Release Bar/Ramp Tension Adjusted

**NOTICE**


When securing the tension on the kevlar cord, make sure the grub screw is on the rivet and NOT on the kevlar cord. The kevlar cord must be beneath the rivet when tightening the grub screw, see [Illustration 15](#) .

Illustration 15: Tightening the Grub Screw

5. Once the Release Bar/Ramp has the appropriate tension applied, tighten the grub screw to secure the tension on the cord. See [Illustration 16](#) .

Illustration 16: Securing the Kevlar Cord

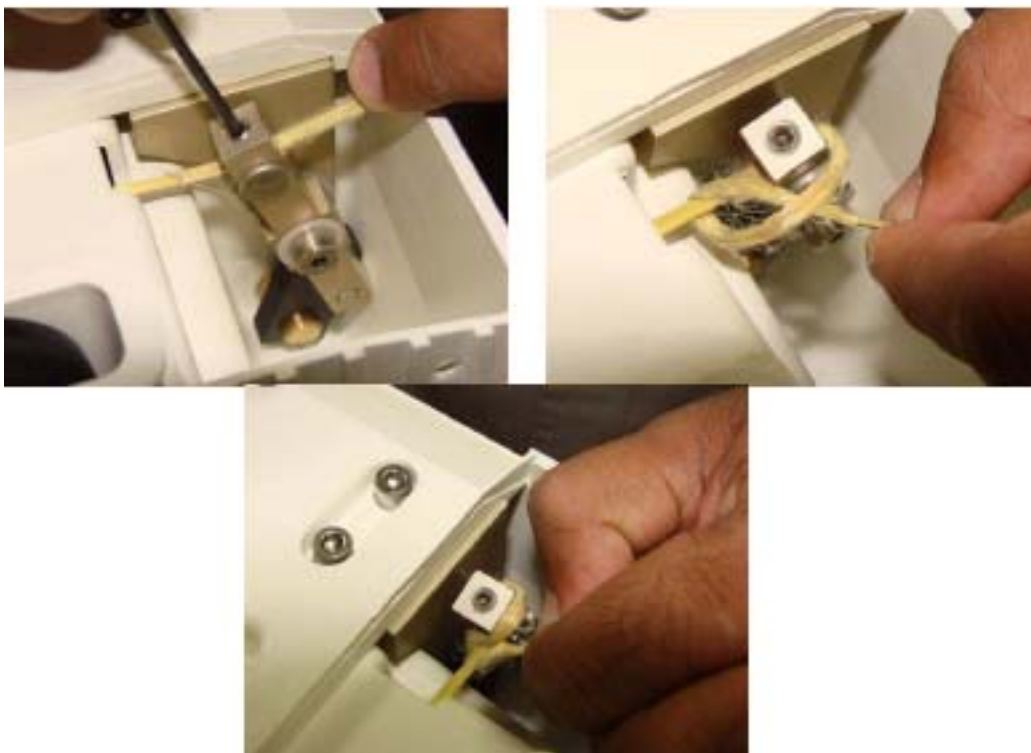6. Tie a knot with the loose end of the kevlar cord as shown above, in [Illustration 16](#) .


7. If the front guide rollers' height needs adjusting, proceed to [Section 1.4.2](#) . If not, then proceed to [Finalization, Step 1](#) .

1.4.2 Front Guide Rollers Height Adjustment

If the Front Guide Rollers require a height adjustment, shims will need to be either added or removed. The Shim part number is 46-317181P1, order additional shims if height is to be raised. The Shims are 0.35 of an inch..

1. If cradle is not removed from the table top, follow the steps in [Step 2](#) . Then proceed to [Section 1.4.1.1](#) through [Section 1.4.1.1, Step 1](#) .
2. Front guide rollers are to be at .25" above the cradle wheel tracks as shown in [Illustration 17](#) . In case of any variation follow the next steps of adjustment process.

Illustration 17: Front Roller Guide Height

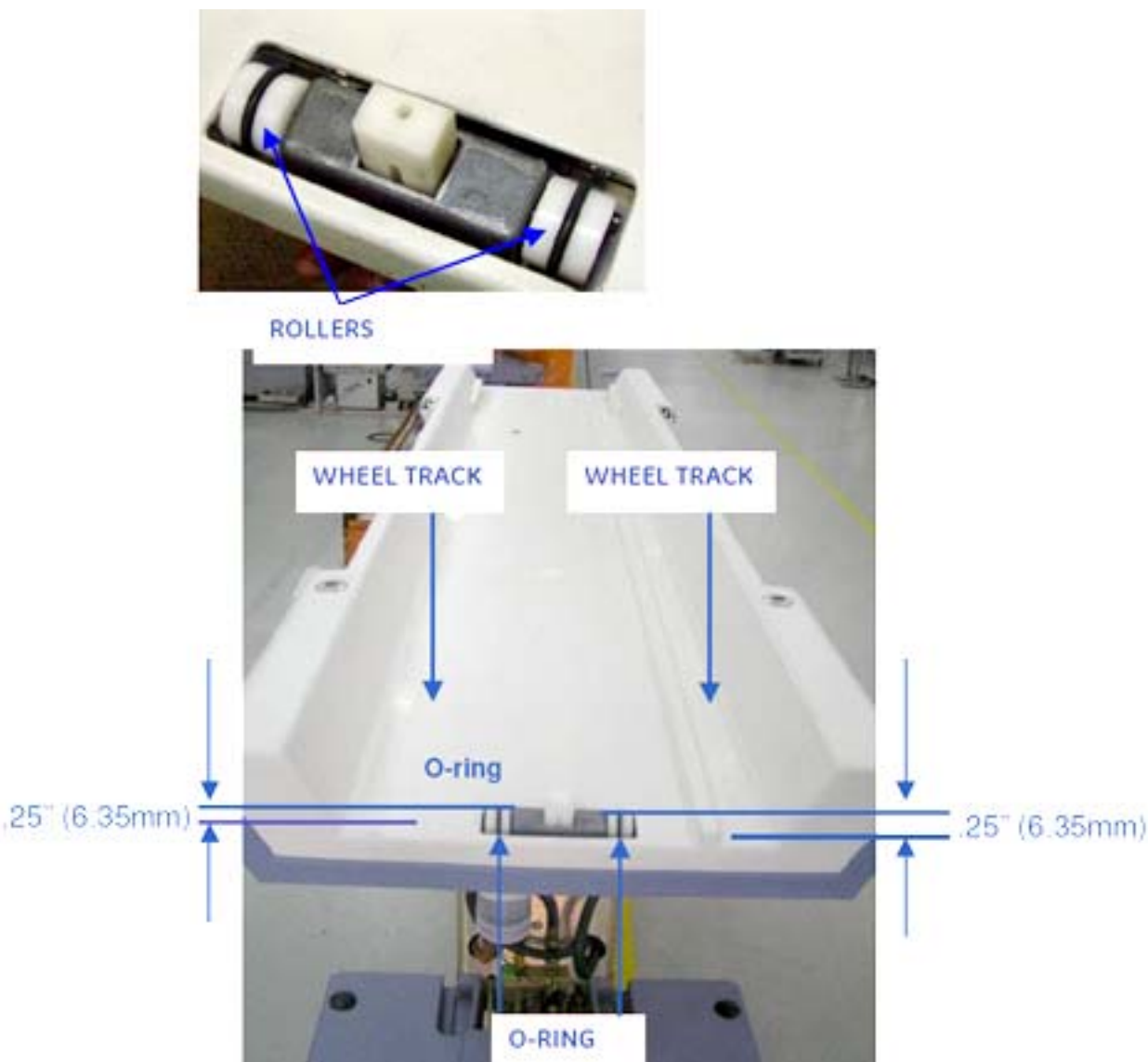3. To readjust the height of the rollers the latch assembly holding blocks must be removed.
See [Illustration 18](#) .

Illustration 18: Location of Holding Blocks

Front latch assembly with guide rollers

Holding block-1

Holding block-2

4. Using the 5/32 Allen Key, loosen the allen screws holding the block, see [Illustration 19](#) .


Illustration 19: Remove the Allen Screws

5. Using the standard screwdriver, remove the screws from the holding block, see [Illustration 20](#) .

Illustration 20: Holding Block Screws

6. Remove the holding blocks from the latch assembly.
7. Add or remove Shims in order to have the O-rings at the correct height of 0.25 inches (6.35 mm) above the horizontal wheel track on either side. See [Illustration 21](#) .

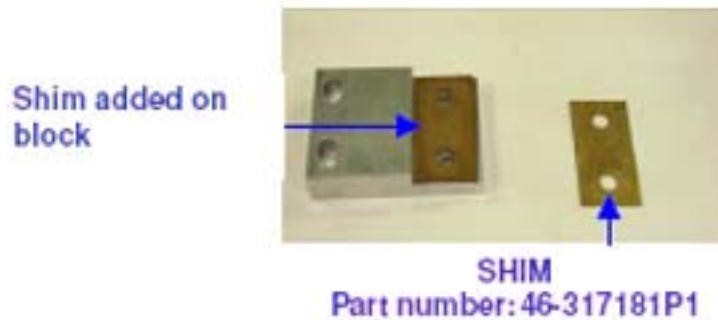Illustration 21: Shimming

8. Once the appropriate height has been achieved, install the holding blocks:
 - a. First, using the standard screwdriver, loosely secure the screws on both holding blocks.
 - b. Secondly, using the 5/32 Allen Key, loosely secure the allen screws on both holding blocks. See [Illustration 22](#) .

Illustration 22: Installing the Holding Blocks

9. Before securing the holding blocks into place, make sure the latch is centered in the slot of the cradle as shown in [Illustration 23](#) .

Illustration 23: Front Latch Position

10. The front latch height needs to be checked before tightly securing the holding blocks in place.
11. Return the Cradle to the table top and put in place to check the latch height.
12. With the casting in the docked position, check that the delrin stopper is flush and the latch is projected 2mm from the casting.
13. Once the latch is at the proper height, secure the screws holding the holding blocks.

1.5 Finalization

1. Re-assemble the cradle by placing the cradle in its position on the table top.
2. Connect the table to the scanner.
3. Move the cradle into the bore of the scanner and then return to its normal position to make sure the cradle is properly placed on the table top.
4. Perform the Emergency Cradle release procedure, by twisting handle to release.
 - a. If it does not release adjust the tension again.
 - b. If the Cradle release functions properly, proceed to [Step 5](#) .

5. Re-engage the cradle to its functioning position.

© 2007 General Electric Company.
GE Medical Systems, a General Electric Company, going to market as GE Healthcare.

www.gehealthcare.com