

MAGNETOM

Installation and Startup

MAGNETOM Option

SW Upgrade VE11

Document Version

Siemens reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order a current version or refer to our website <http://www.healthcare.siemens.com>.

Disclaimer

Siemens provides this documentation “as is” without the assumption of any liability under any theory of law.

The content described herein requires superior understanding of our equipment and may only be performed by qualified personnel who are specially trained for such installation and/or service.

Copyright

“© Siemens, 2014” refers to the copyright of a Siemens entity such as:

Siemens Healthcare GmbH - Germany

Siemens Aktiengesellschaft - Germany

Siemens Shenzhen Magnetic Resonance Ltd. - China

Siemens Shanghai Medical Equipment Ltd. - China

Siemens Healthcare Private Ltd. - India

Siemens Medical Solutions USA Inc. - USA

Siemens Healthcare Diagnostics Inc. - USA and/or

Siemens Healthcare Diagnostics Products GmbH - Germany

1	General Information	4
1.1	General Information	4
1.1.1	Purpose	4
1.1.2	New SW Features	4
1.1.3	syngo Expert-i	4
1.1.4	Prerequisites	5
1.1.5	Special Tools/Documents	6
1.2	Material Information	7
1.2.1	Ordering Information	7
1.2.2	Content of the Update Kit	7
1.2.3	Return of Parts	8
2	CSSI Upgrade First Steps for VE11	9
2.1	CSSI Upgrade First Steps for VE11A/VE11B	9
2.1.1	Work Steps for creating the backup	9
2.2	CSSI Upgrade First Steps for VE11C	32
2.2.1	Work Steps for creating the backup	32
3	CSSI Upgrade Second Steps for VE11	57
3.1	CSSI Upgrade Second Steps for VE11A/VE11B	57
3.1.1	Work Steps	57
3.2	CSSI Upgrade Second Steps for VE11C	114
3.2.1	Work Steps	115
4	Additional Information	178
4.1	Appendix	178
4.1.1	Drive Letter Assignment for the CD/DVD ROM Drive and/or for the CD/DVD-R Drive	178
4.1.2	Correct syngoUI Language Settings	180
4.1.3	Failure messages where the regional language setting did not match, e.g. UI language	182
4.1.4	Language settings	184
4.1.5	Deactivation of OLM	216
4.1.6	New Password Handling in syngo	227
4.1.7	Advanced User handling	231
4.1.8	New feature for DICOM node configuration	233
4.1.9	System management installation workaround	234
4.1.10	Memory upgrade to 32 GByte	241
5	Changes to previous Version	249
6	List of Hazard IDs	250

1.1 General Information

1.1.1 Purpose

This SW Update provides the current software syngo MR E11

1.1.2 New SW Features

The syngo MR E11 provides the new Dot Cockpit for easy Dot Engine configuration together with several workflow and performance enhancements. There are new licensed options available with the syngo MR E11 software:

- Quiet Suite
- MyoMaps
- TWIST VIBE
- StarVIBE
- Advanced WARP
- LiverLab

1.1.3 syngo Expert-i

syngo Expert-i

- Interactive remote assistance during the MR exam. Get real-time access to imaging data and exam information from any networked PC within the hospital network. Until now, radiologists or other experts had to stop what they were doing and go to the MR scanner when complex clinical questions arose.

Now, questions can be addressed quickly and efficiently via a networked PC.

Some of the benefits of syngo Expert-i are:

- » Excellent results right from the first examination
- » Streamlined workflow & faster patient throughput
- » Reduces repeat rates with a check on images while the patient is still in the examination room
- » Reduces training effort by enabling expert assistance for specialized procedures

1.1.4 Prerequisites



Every CSSI upgrade package with a VE11A COBR-Recovery-Image DVD creation date of 24/07/2015 and newer includes the OLM deactivation.

In case one of the following UI's:

- MR019/15/R with VB13A-SP07
- UI MR021/15/R with VD13C-SP03

was not performed, the OLM is still active. You have to apply the OLM deactivating description according Appendix (→ Deactivation of OLM / Page 216)

In case the OLM is already deactivated, one of the mentioned UI's was performed. No further actions are necessary.

In case your VE11A COBR-Recovery-Image DVD presents a creation date before of 24th July 2015 and the OLM was already deactivated you have to perform this UI together with UI MR005/15/P.

Refer also to the Service Knowledge Base entry: SKB0067279.



"Work In Progress" packages provided for a previous software version (VD13A/VD13B/VD13C/VD13D/VD13E) will generally NOT work with VE11A or VE11B or VE11C!



The product brochure **DotGO Quick Guide** (MR-05022.640.01.02.xx) as part of the ADOC has to be handed over to the customer. It includes useful information about the new software.



Due to DICOM changes associated with the VE11A/VE11B/VE11C, adjustments involving the locally connected DICOM nodes may become necessary on site.



The BreVis evaluation SW is already included in the VE11A/VE11B/VE11C. The ordering of a separate installation medium is no longer necessary.

- Before the VE11A SW installation the UI MR016/14/P has to be performed on the Skyra systems with the serial number 45179 to 45590 with the exception: 45378, 45397 and 452424.
- To also use IDEA with MR systems running VE11A/VE11B, a current version of the IDEA software is necessary.
 - Request the current IDEA package from H IM MR PI CM for sites where this feature is being used.
 - Please note: In accordance with the Siemens policy, IDEA will be provided within the framework of a cooperation agreement only. In case of questions, please contact MR cooperation management - H IM MR PI CM.

- In case a Leonardo/MMWP workstation (MultiModality Workplace) is used to process data acquired with the MR system with NUMARIS/4 E11 the current Leonardo/MMWP workstation (MultiModality Workplace) versions are not compatible. Required is the new Leonardo/MMWP workstation (MultiModality Workplace) version MM WP VE61A. The release of this version is scheduled for March 2015.

The following host PC types are supported by the MMWP version VE61A:

- R670-2
- M470-2
- M720

If the customer has a different host PC the hardware must be exchanged too.

- Ensure that the VE11A/VE11B/VE11C application training timing is suitable for your site.
- One of the following UI's for disable the OLM must be performed:
 - for MR systems running VD13A: UI MR019/15/R
 - for MR systems running VD13B: UI MR020/15/R
 - for MR systems running VD13C: UI MR021/15/P
 - for MR systems running VD13D: UI MR022/15/R
 - for MR systems running VD13E: UI MR023/15/R



Application Training required!

1.1.5 Special Tools/Documents

Please check for the latest available version of this document in the Intranet or contact your local Siemens office if you do not have access to the Siemens Intranet/Extranet.

- Anti-static wristband and grounded tabletop pad
(TDK equipment, Type 8005 from 3M Deutschland GmbH)

1.2 Material Information

1.2.1 Ordering Information


The VE11A/VE11B SW material has to be ordered via the PM process.

Quantity	Update kit
1	14441798

1.2.2 Content of the Update Kit

Material No.	Name	MRAWP only	MRAWP + MRWP
10928179	BIOS CD syngo MR VE11-Container	1	1
10432853	DVD COBR-Recovery-Image	1	2
10915873	Expert-I	1	1
07752103	Set of empty DVDs	1	2
	Customer documentation package	1	1
	Installation & Start-up: SW Upgrade VE11A/VE11B (Option: MR00-000.814.56.xx.02)	1	1
	Customer Information Kit (variabel)	0 or 1	0 or 1
10928151	MEMORY MODULE 4GB	4	4
11234279	Retrofit measuring system	n.a.	n.a.

Fig. 1: Content of the retrofit measuring system kit, p. n. 11234279

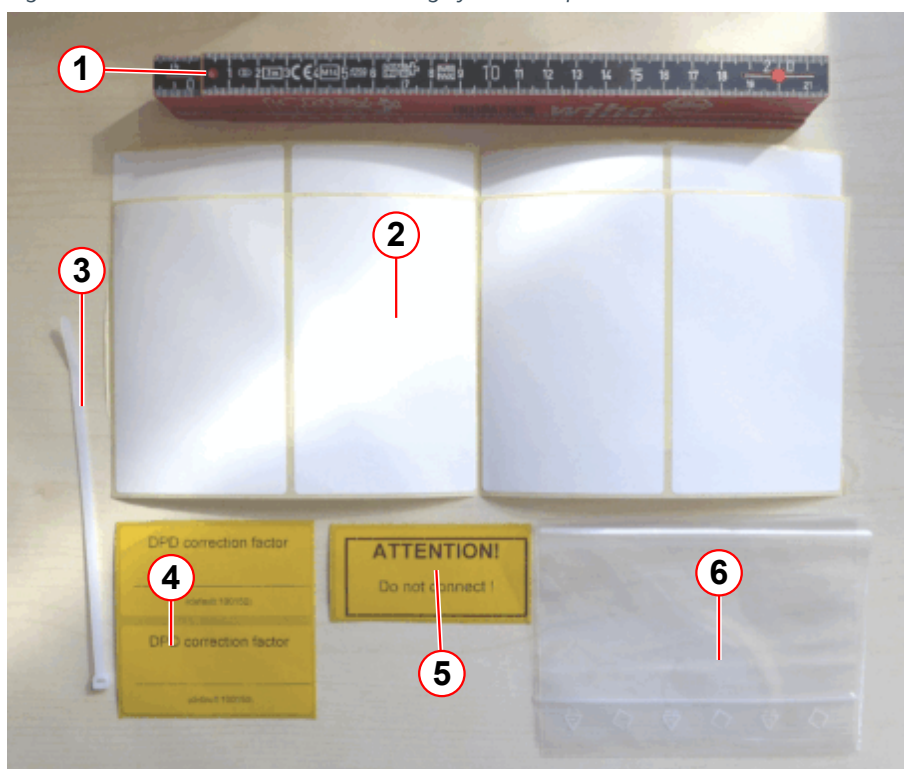

Pos.	Quantity	Material No.	Name
1	1	10759380	Folding ruler, 2 m (EU standard accuracy class 2)
2	8	n.a.	Set of white tape strips, 7.4 x 10.5 cm
3	1	n.a.	Tie wrap, 18.8 X 0.48 cm
4	2	n.a.	"DPD correction factor" label, 7.5 x 3.5 cm
5	1	n.a.	"Attention! Do not connect" label, 7.5 x 3.5 cm
6	1	n.a.	Plastic bag

1.2.3 Return of Parts

n.a.

2.1 CSSI Upgrade First Steps for VE11A/VE11B



The following steps describes the procedure which creates the backup for the upgrade to VE11A/VE11B. This backup has to be send for the CSSI creation.

2.1.1 Work Steps for creating the backup



If you have a satellite console with a shared database connected to the main console, perform the following steps first.

If no MRSC is present, skip the following steps and start at (→ MRAWP Step: Backing up dynamic Data at the MRAWP (MRC) / Page 12)


2.1.1.1 Overview

Performing the Backup for the CSSI Update on-site




Please start the backup of the main console MRC directly after the start of the backup creation in the local service UI.

Fig. 2: Flowchart backup on-site

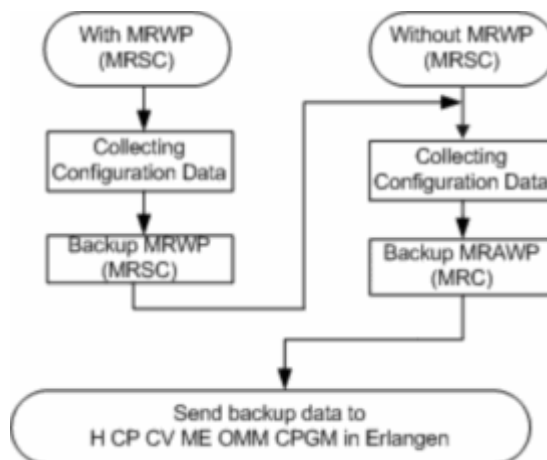

Performing a Remote Backup for the CSSI Update

Fig. 3: Flowchart - remote backup

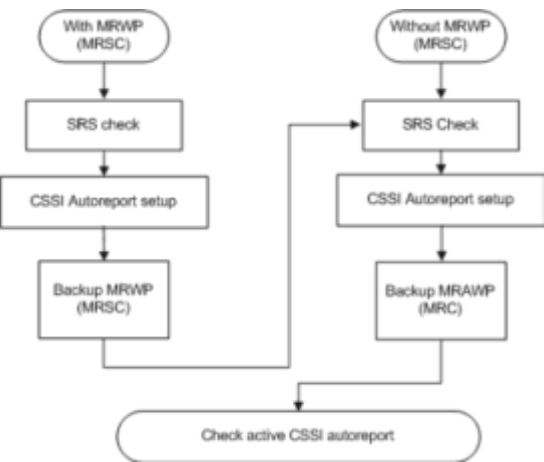

2.1.1.2 Perform on-site backup for the CSSI update

MRWP Step: Backing up Dynamic Data at the MRWP (MRSC)

Collecting MRWP Configuration Data from Local Service

ACTION	MRAWP
■ Select Options --> Service --> Local Service.	☐
■ Enter the valid service key and click OK.	☐
■ Select Reports in the service platform.	☐
■ On the Report Files page, select Site Info Imprint	☐
■ Print out the Imprint report (generate a PDF file or write down the system options) Launch the update of the Imprint file (bottom right corner) In order to generate a PDF file of the Imprint report, press " show as PDF ". Save the PDF file under C:\MedCom\MriSiteData\Imprint.pdf	☐
■ In the Report Files page, select Site Info Configuration	☐

ACTION	MRAWP
- Print out the Configuration report (generate a PDF file or write down the system options) - DICOM networking (write down all DICOM network nodes and settings) - DICOM print devices (write down all DICOM print devices and settings) - DICOM HIS/RIS (write down all DICOM HIS/RIS nodes and settings) In order to generate a PDF file of the Configuration report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Report.pdf	☐
In the Report Files page, select Site Info MR Configuration	☐
Print out the MR Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MrReport.pdf	☐
Click EXIT to close the Report Files popup window	☐



It is recommended to copy the PDF Report Files to an external medium to have the information for the configuration check in the new software version.

Please check the external media first with the virus scanner of the service PC. Be sure that the external medium used is free of viruses.

Masterbackup at the MRWP (MRSC)

ACTION	MRSC
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Systemserial# Master Backup VD13A or VD13B or VD13C or VD13D or VD13E MRWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK.	☐

ACTION	MRSC
Note: ■ In some cases, the following popup message: "Please insert the media into Drive S: \ and press so that 'OK' to proceed appears. Click "OK" to continue ■ If the following error messages are displayed at the end of the backup routine, ignore them. The backup is then completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
■ Remove the CD.	☐
■ Click Home and close the local service UI.	☐



If the backup creation to CD fails, please create the master backup to harddisc and transfer it to your service PC. On your service PC burn the master backup to CD.

MRAWP Step: Backing up dynamic Data at the MRAWP (MRC)






If you have a satellite console with a shared database connected to the main console, start with the action item:

- (→ MRWP Step: Backing up Dynamic Data at the MRWP (MRSC) / Page 10)

This item has to be performed first. If this step has already been completed on the MRSC, continue with the MRC.

Collecting MRAWP Configuration Data from Local Service

ACTION	MRWP
■ Select Options --> Service --> Local Service.	☐
■ Enter the valid service key and click OK.	☐
■ Select Reports in the service platform.	☐

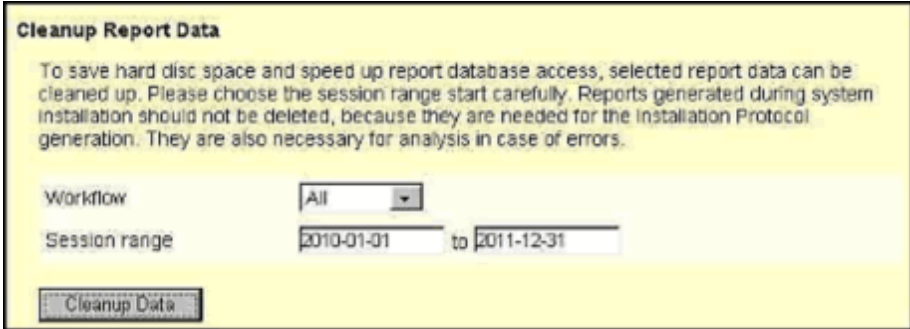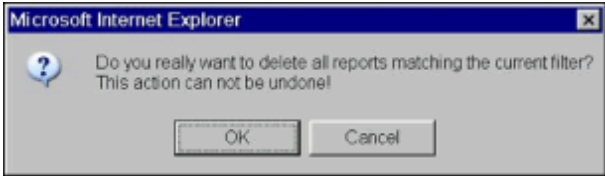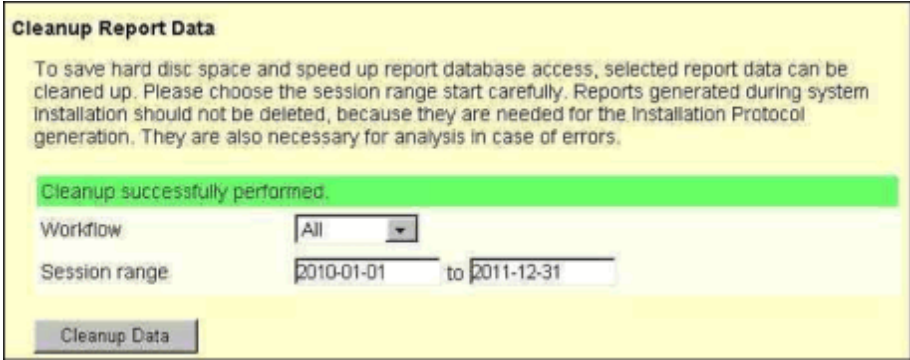ACTION	MRWP
■ Purge service-related files (test tool files and tune-up, QA files older than one year). If the files are not purged, the size of the created Master backup will be huge. - Select Reports ->Clean up - Under Cleanup Report Data select Workflow All - Session range 2010-01-01 to (today - 1 year) e.g. 2011-07-01 - Also you can click in the fields and a calendar popup appears. <i>Fig. 4: Data Cleanup</i>  (1) Start date (2) End date - Press the Cleanup Data button. <i>Fig. 5: Confirmation popup</i>  - Confirm the popup with OK. - The following shows a successful file purge. <i>Fig. 6: Successful cleanup message</i> 	☐
■ On the Report Files page, select Site Info Imprint	☐


ACTION	MRWP
Print out the Imprint report (generate a PDF file or write down the system options) Launch the update of the Imprint file (bottom right corner) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Imprint.pdf	☐
In the Report Files page, select Site Info Configuration	☐
- Print out the Configuration report (generate a PDF file or write down the system options) - DICOM networking (write down all DICOM network nodes and settings) - DICOM print devices (write down all DICOM print devices and settings) - DICOM HIS/RIS (write down all DICOM HIS/RIS nodes and settings) In order to generate a PDF file of the Configuration report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Report.pdf	☐
In the Report Files page, select Site Info MR Configuration	☐
Print out the MR Configuration Report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MrReport.pdf.	☐
In the Report Files page, select Site Info CAN24 Configuration	☐
Print out the CAN24 Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\CAN24Report.pdf	☐
In the Report Files page, select Site Info MBox Configuration	☐
Print out the MBox Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MBoxReport.pdf	☐
Click EXIT to close the Report Files popup window	☐



It is recommended to copy the PDF report files to an external medium to have the information for the configuration check in the new software version.

Please check the external media first with the virus scanner of the service PC. Be sure that the external medium used is free of viruses.

Masterbackup at the MRAWP (MRC)

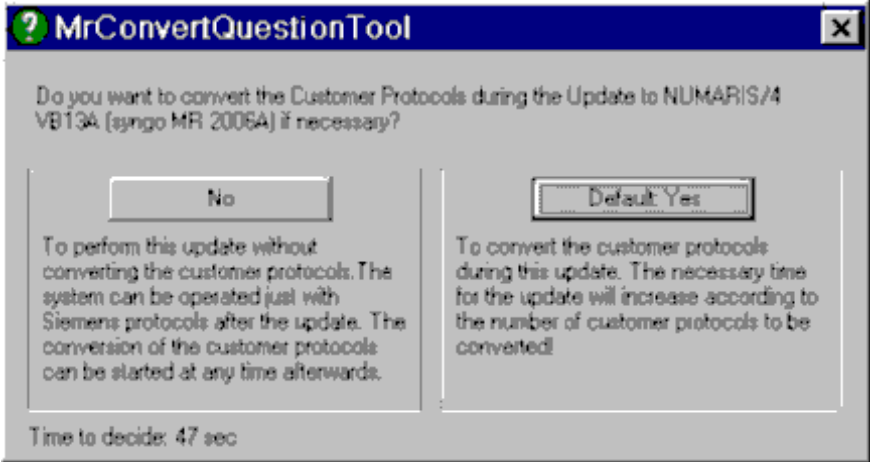ACTION	MRAWP
■ Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Master Backup VD13A or VD13B or VD13C or VD13D or VD13E MRAWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK. During the Master Backup, the following popup appears: Fig. 7:  Press the "Yes" button. At this point, this popup will not start the protocol conversion if "yes" was selected. This popup determines the default focus ("Yes" or "NO") of the MrConvertQuestion-Tool when restoring the Master Backup.	<input data-bbox="1396 582 1428 622" type="checkbox"/>

ACTION	MRAWP
Note: ■ In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed," appears. Click "OK" to continue ■ If the following error messages are displayed at the end of the backup routine, ignore them. The backup is then completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed .. and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ..	
■ Remove the CD.	☐
■ Click Home and close the local service UI.	☐



If the backup creation to CD fails, please create the master backup to harddisc and transfer it to your service PC. On your service PC burn the master backup to CD.

Sending the Backup CDs to CV ME

Send the backup CDs to the following address:

Siemens Healthcare GmbH

HC CP CV ME OMM CPGMC

Attn.: Mr. G. Pfaendner

Allee am Roethelheimpark 2

D - 91052 Erlangen

"SW Sales / CSSI Customer Data"

2.1.1.3 Perform remote backup for the CSSI update



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the request for full access. Also inform the customer to keep the system running in standby mode for the backup transfer.

Prerequisites

- System is connected to the SRS infrastructure via a VPN connection

- System is in **full remote access mode**
- The estimated time is about 30-40 min for the remote backup, depending on the amount of backup data.

SRS Configuration Check

- Remote login to the system using the cRSP (common Remote Service Platform)
- Check the RDIAG function on the **MRAWP** and, if applicable, on the **MRWP** as well
 - Test the FTP connection to the RDIAG server as follows:
 Local Service -> File&Image Tools
 Source: File
 Action: Transfer to System
 File Type: Any file
 From source: SRS_ftp
 If a list from the directory /outgoing/mr (directory on the RDIAG server) is displayed, the connection and FTP access are OK



If the FTP test fails, check the SRS configuration; refer to (→ Final Check - to be performed for MRC and MRSC (if applicable) / Page 110)



If you do not have a satellite console with a shared database connected to the main console, the step (→ MRWP Step: Auto Report setup and remote backup / Page 17) can be skipped.

Continue then with (→ MRAWP Step: Auto Report Setup and remote Backup / Page 23).

MRWP Step: Auto Report setup and remote backup



To create the backup, full access to the system is required. Please inform the customer to get full access. The time required for full access mode is about 30 minutes.

Setup of the CSSI Auto Report at the MRWP




The setup of the CSSI autoreport is only necessary, if they are not available on the system. Please check if the CSSI autoreports are enabled.

ACTION	MRWP
■ Login to the system via cRSP (common Remote Service Platform).	☐
■ Enter the valid service key and click OK .	☐
■ Select Auto Report in the service platform.	☐



ACTION	MRWP
Under Report, select the entry <code>TEMPLATE</code>.	☐
- Create the new Auto Report with the name <code>CSSI_Tran_N4</code> Use the following values: - Active: Mark it - Name: <code>CSSI_Tran_N4</code> - Target: <code>SRS_mail</code> - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: <code>C:\MedCom\service\transfer*Numaris4*.ar</code> - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: <code>SRS_ftp</code> - Send when: At once	☐




Fig. 8: `CSSI_Tran_N4`

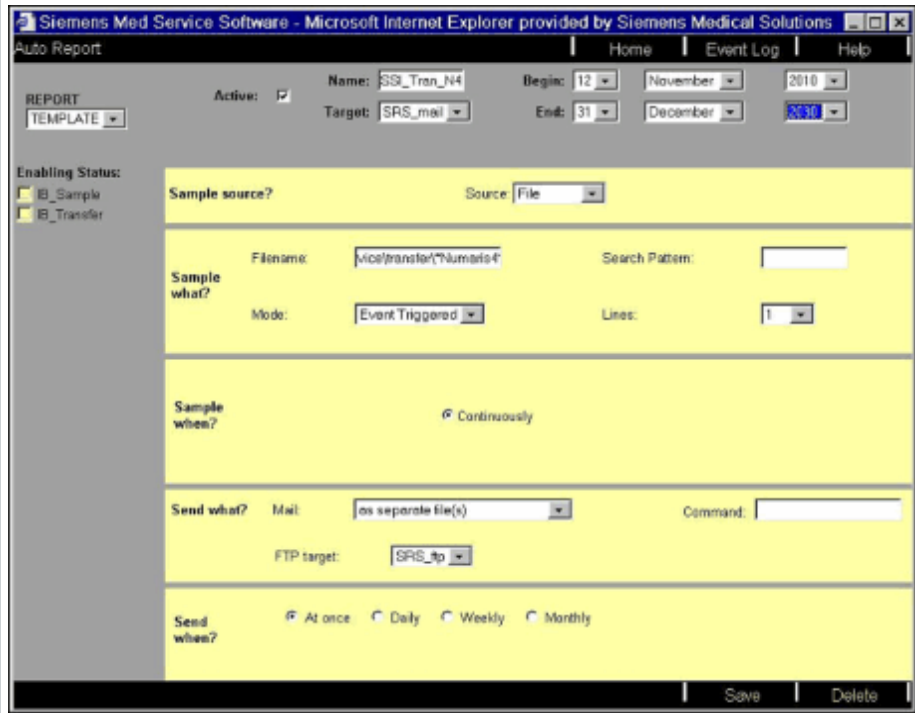

- Click **Save**.

☐

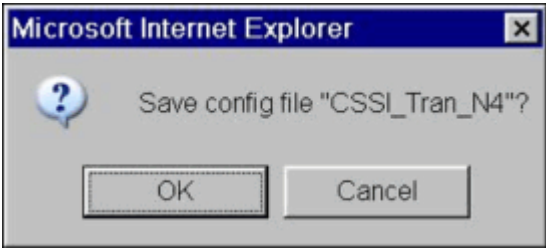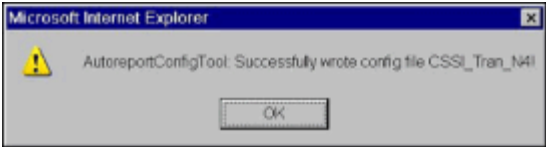ACTION	MRWP
Confirm the pop-up with OK. Fig. 9: Save OK 	☐
Continue by clicking OK on the pop-up. Fig. 10: Success pop-up 	☐
Under Report, select the entry TEMPLATE.	☐
- Create the new Auto Report with the name CSSI_Tran_SecSe Use the following values: - Active: Mark it - Name: CSSI_Tran_SecSe - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*Security-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

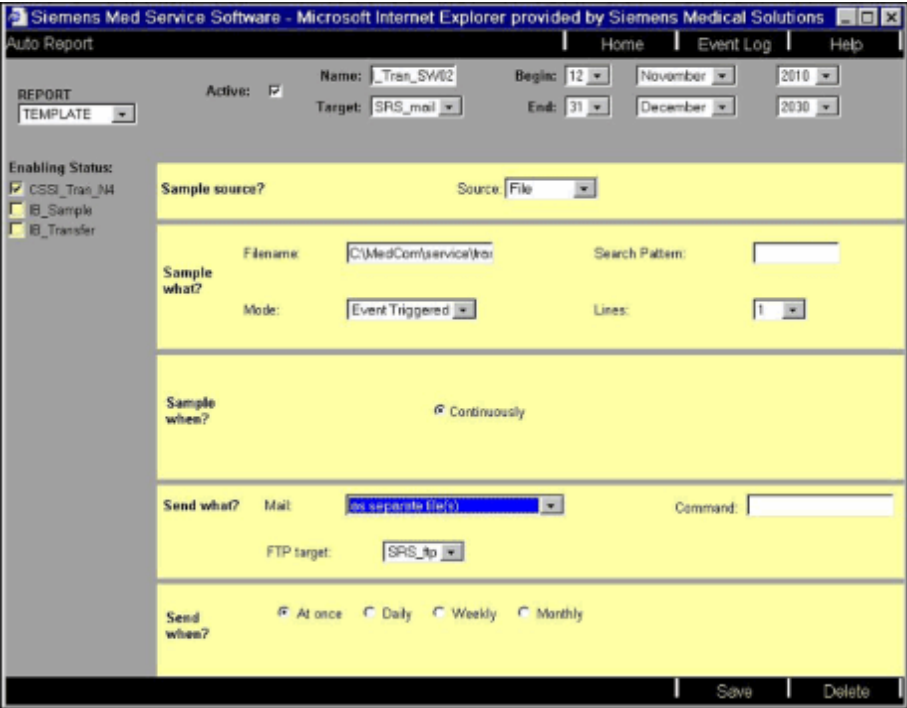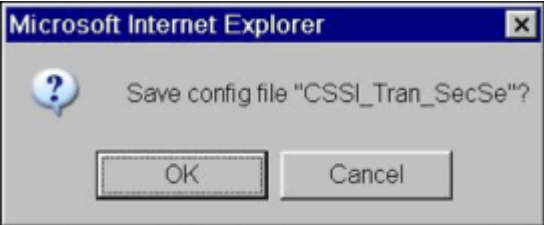
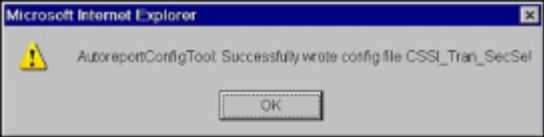
ACTION		MRWP
Fig. 11: Autoreport CSSI_Tran_SecSe 		
■ Click Save .		☐
■ Confirm the pop-up with OK .		☐
Fig. 12: Save OK 		
■ Continue by clicking OK on the pop-up.		☐
Fig. 13: Success pop-up 		
■ Under Report , select the entry TEMPLATE .		☐

ACTION	MRWP
■ Create the new Auto Report with the name CSSI_Trans_SW02 Use the following values: - Active: Mark it - Name: CSSI_Trans_SW02 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*SW-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

Fig. 14: Autoreport CSSI_Trans_SW02

- Click **Save**.

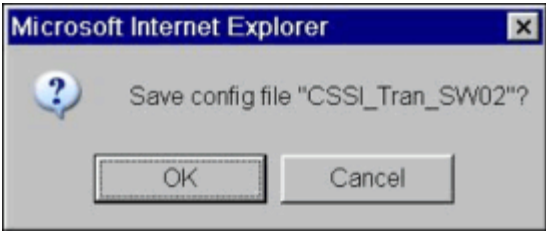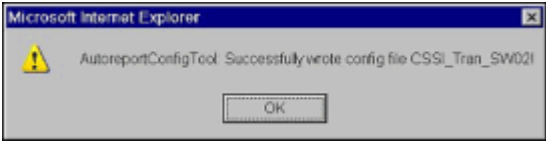

ACTION	MRWP
Confirm the pop-up with OK. <i>Fig. 15: Save OK</i> 	☐
Continue by clicking OK on the pop-up. <i>Fig. 16: Success pop-up</i> 	☐
Check the newly created Autoreports to see if they are configured correctly and are activated.	☐
Click Home to leave the Auto Report section.	☐

Performing the Backup at the MRWP



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the request for full access. Also inform the customer to keep the system running in standby mode for the backup transfer.

ACTION	MRWP
Select Backup&Restore in the service platform.	☐
- Perform the backup in local service. The backup has to be performed for the following packages: Numaris4, SW settings02, and Security Settings. Select Backup in the Command menu, select [-E-] FIXED in the Drives menu. - Select SW Settings02 in the Package menu, start the backup by clicking OK. - Select Security Settings in the Package menu, start the backup by clicking OK. - Select Numaris4 in the Package menu, start the backup by clicking OK. Check for errors using the log window.	☐
Click Home to leave Backup & Restore	☐

Move Backups into the transfer directory MRWP

ACTION	MRWP
Select Utilities in the service platform.	☐
Under Source, select the entry Escape to OS in the pull-down menu	☐
Under Command, keep the default entry NT Command Interpreter.	☐
Under Parameters, type in the following: <pre>move /y E:*.ar C:\MedCom\service\transfer\</pre>	☐
Click Go. The three backups that have been created are now moved to the directory C:\Med-Com\service\transfer.	☐
Fig. 17: Move output 	
Click Home to leave Utilities	☐
Check the Auto Reports to see whether they are active or not (see (→ Checking active CSSI Auto Reports / Page 29))	☐
Close the service UI in order to close the connection.	☐

MRAWP Step: Auto Report Setup and remote Backup

Setup of the CSSI Auto Report at the MRAWP





The setup of the CSSI autoreport is only necessary, if they are not available on the system. Please check if the CSSI autoreports are enabled.




ACTION	MRAWP
Login to the system via cRSP (common Remote Service Platform).	☐
Select Auto Report in the service platform.	☐
Under Report, select the entry TEMPLATE.	☐
- Create the new Auto Report with the name CSSI_Trans_N4 Use the following values: - Active: Mark it - Name: CSSI_Trans_N4 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*Numaris4*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

Fig. 18: CSSI_Trans_N4

- Click **Save**.

☐

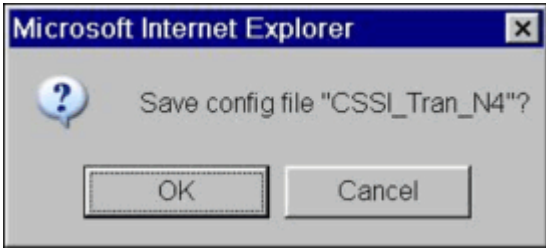
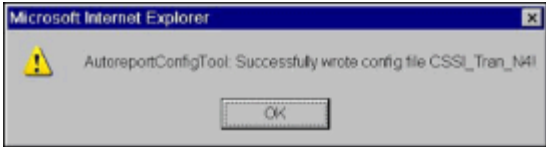
ACTION	MRAWP
Confirm the pop-up with OK. Fig. 19: Save OK 	☐
Continue by clicking OK on the pop-up. Fig. 20: Success pop-up 	☐
Under Report, select the entry <code>TEMPLATE</code>.	☐
- Create the new Auto Report with the name <code>CSSI_Tran_SecSe</code> Use the following values: - Active: Mark it - Name: <code>CSSI_Tran_SecSe</code> - Target: <code>SRS_mail</code> - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: <code>C:\MedCom\service\transfer*Security-Settings*.ar</code> - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: <code>SRS_ftp</code> - Send when: At once	☐

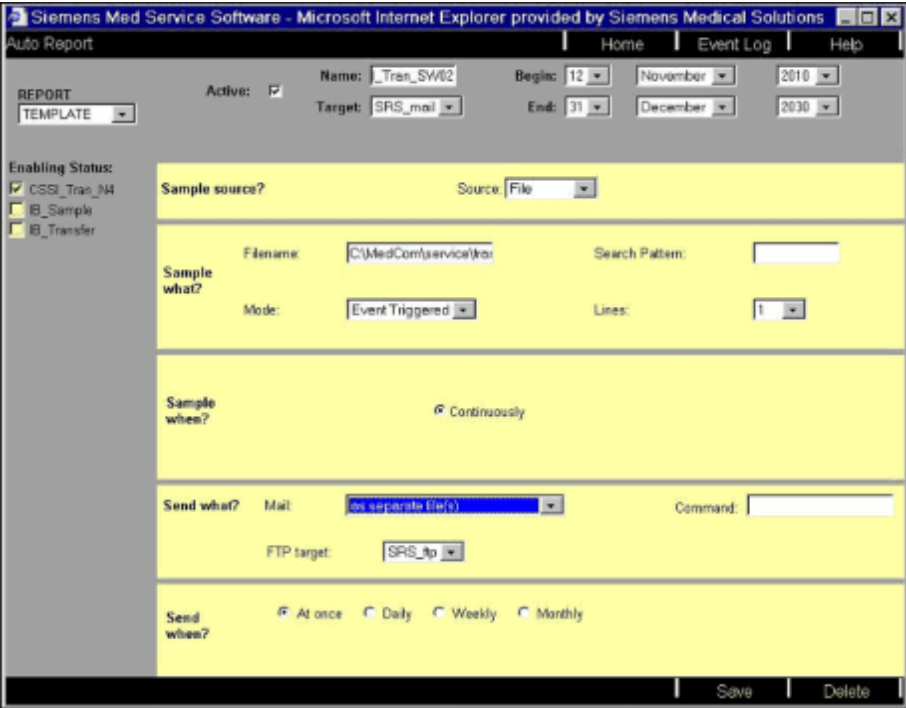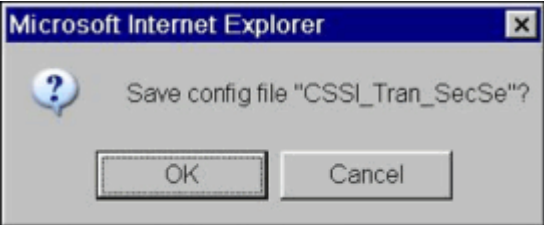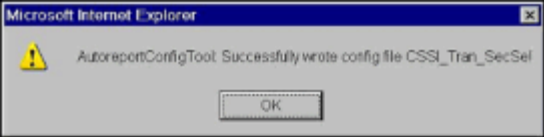ACTION		MRAWP
Fig. 21: Autoreport CSSI_Tran_SecSe 		
■ Click Save .		☐
■ Confirm the pop-up with OK .		☐
Fig. 22: Save OK 		
■ Continue by clicking OK on the pop-up.		☐
Fig. 23: Success pop-up 		
■ Under Report , select the entry TEMPLATE .		☐

ACTION	MRAWP
■ Create the new Auto Report with the name CSSI_Trans_SW02 Use the following values: - Active: Mark it - Name: CSSI_Trans_SW02 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*SW-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

Fig. 24: Autoreport CSSI_Trans_SW02

■ Click Save.	☐
--	--------------------------

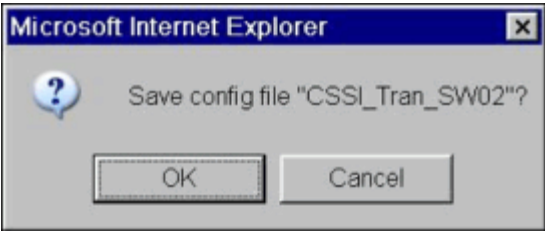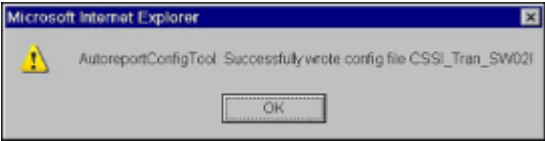ACTION	MRAWP
Confirm the pop-up with OK. Fig. 25: Save OK 	☐
Continue by clicking OK on the pop-up. Fig. 26: Success pop-up 	☐
Check the newly created Autoreports to see if they are configured correctly and are activated.	☐
Click Home to leave the Auto Report section.	☐

Performing the Backup at the MRAWP



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the request for full access. Also inform the customer to keep the system running in standby mode for the backup transfer.



ACTION	MRAWP
Select Backup&Restore in the service platform.	☐
- Perform the backup in local service. The backup has to be performed for the following packages: Numaris4, SW settings02, and Security Settings. Select Backup in the Command menu, select [-E-] FIXED in the Drives menu. - Select SW Settings02 in the Package menu, start the backup by clicking OK. - Select Security Settings in the Package menu, start the backup by clicking OK. - Select Numaris4 in the Package menu, start the backup by clicking OK. The MRConvertQuestion pop-up is not displayed when backup is performed via remote login. The duration of the Numaris4 backup is about 15-20 min depending on the amount of data. Check for errors using the log window.	☐
Click Home to leave Backup & Restore	☐


Move Backups into the transfer directory

ACTION	MRC
Select Utilities in the service platform.	☐
Under Source, select the entry Escape to OS in the pull-down menu	☐
Under Command, keep the default entry NT Command Interpreter.	☐
Under Parameters, type in the following: <pre>move /y E:*.ar C:\MedCom\service\transfer\</pre>	☐
Click Go. The three backups that have been created are now moved to the directory C:\Med-Com\service\transfer. Fig. 27: Move output 	☐
Click Home to leave Utilities	☐

Checking active CSSI Auto Reports

ACTION	MRAWP/MRWP
Select File & Image Tools in the service platform.	☐
Under Source, select the entry Transfer List from the pull-down menu.	☐
Under Action, select the entry Show/Delete in the pull-down menu.	☐
Now the Transfer List is shown. Check to see if the FTP transfer for the backup files has been started. For that purpose, select the FTP entries in the transfer manager list.	☐

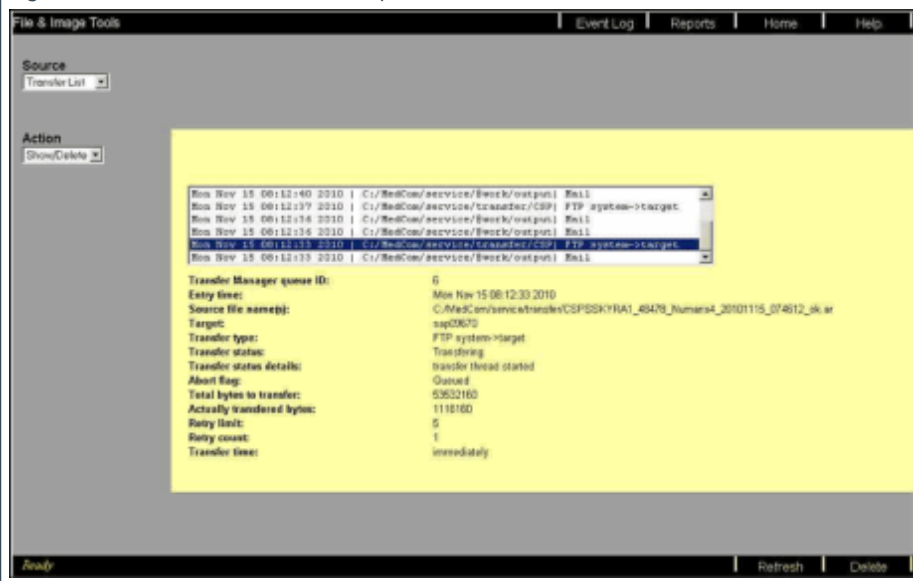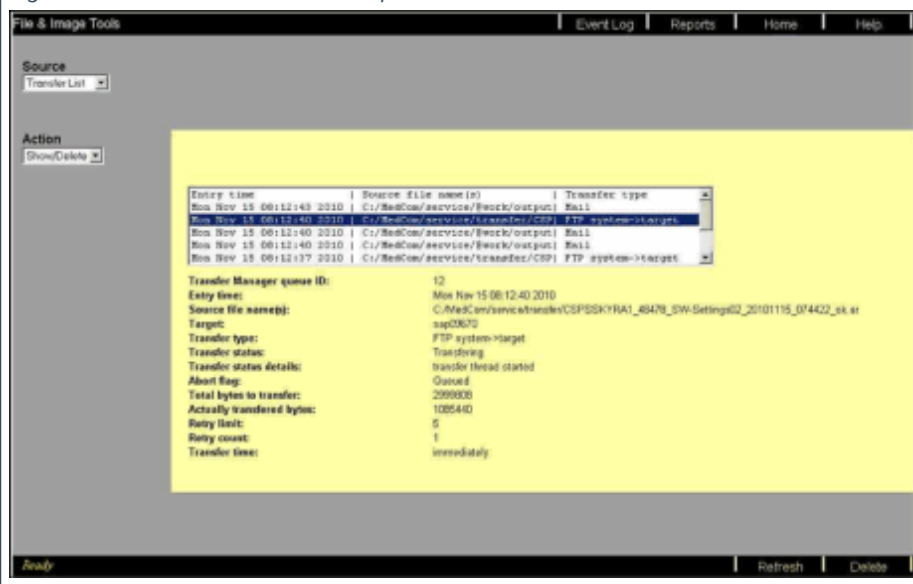ACTION	MRAWP/MRWP
■ For the Numaris4 backup file. Fig. 28: Transfer List active autoreport Numaris4 	☐
■ For the Software02 backup file. Fig. 29: Transfer List active autoreport Software02 	☐

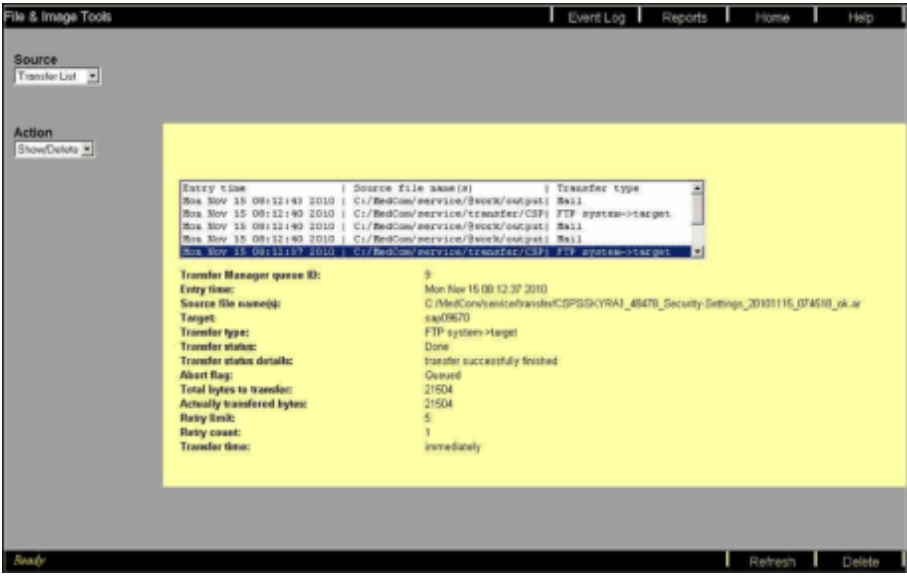ACTION	MRAWP/MRWP
For the Security Settings backup file. <i>Fig. 30: Transfer List active autoreport SecuritySettings</i> 	☐
Click Home to leave Utilities	☐
Close the local service UI.	☐



Please inform the customer not to switch off the overall system. The system should be switched to standby.

2.2 CSSI Upgrade First Steps for VE11C



The following steps describe the procedure which creates the backup for the upgrade to VE11C. This backup has to be sent for the CSSI creation.

2.2.1 Work Steps for creating the backup



If you have a satellite console with a shared database connected to the main console, perform the following steps first.

If no MRSC is present, skip the following steps and start at (→ MRAWP Step: Backing up Dynamic Data at the MRAWP (MRC) / Page 35)


2.2.1.1 Overview

Performing the Backup for the CSSI Update on-site




Please start the backup of the main console MRAWP (MRC) directly after starting the backup creation in the local service UI.

Fig. 31: Flowchart backup on-site

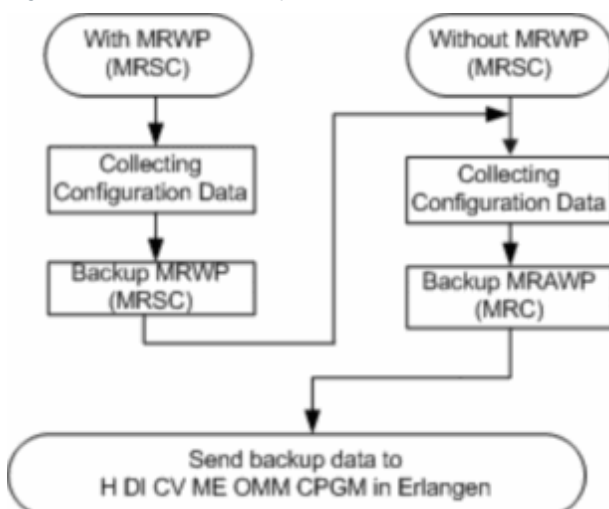

Performing a Remote Backup for the CSSI Update

Fig. 32: Flowchart - remote backup

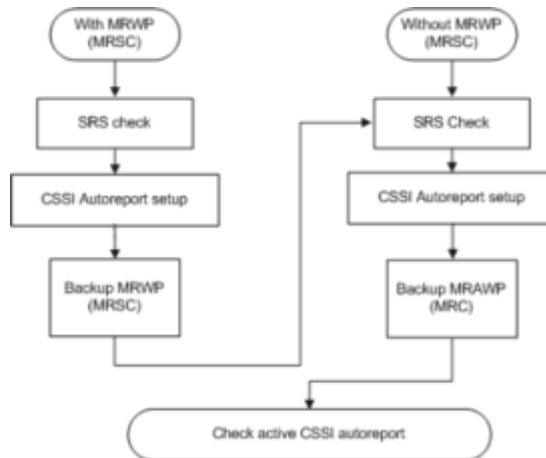

2.2.1.2 Perform on-site backup for the CSSI update

MRWP Step: Backing up Dynamic Data at the MRWP (MRSC)

Collecting MRWP Configuration Data from Local Service

ACTION	MRAWP
■ Select Options --> Service --> Local Service .	☐
■ Enter the valid service key and click OK .	☐
■ Select Reports in the service platform.	☐
■ On the Report Files page, select Site Info Imprint	☐
■ Print out the Imprint report (generate a PDF file or write down the system options) Launch the update of the Imprint file (bottom right corner) In order to generate a PDF file of the Imprint report, press " show as PDF ". Save the PDF file under C:\MedCom\MriSiteData\Imprint.pdf	☐
■ On the Report Files page, select Site Info Configuration	☐

ACTION	MRAWP
- Print out the Configuration report (generate a PDF file or write down the system options) - DICOM networking (write down all DICOM network nodes and settings) - DICOM print devices (write down all DICOM print devices and settings) - DICOM HIS/RIS (write down all DICOM HIS/RIS nodes and settings) In order to generate a PDF file of the Configuration report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Report.pdf	☐
On the Report Files page, select Site Info MR Configuration	☐
Print out the MR Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MrReport.pdf	☐
Click EXIT to close the Report Files popup window	☐



It is recommended to copy the PDF Report Files to an external medium to have the information for the configuration check in the new software version.

Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

Master Backup at the MRWP (MRSC)

ACTION	MRWP
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Systemserial# Master Backup VD13 or VE11A or VE11B MRWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK.	☐

ACTION	MRWP
Note: ■ Check the output log for errors. If any other error occurs apart from those mentioned below, the master backup has to be created again. ■ In some cases, the following popup message: "Please insert the media into Drive S: \ and press 'OK' to proceed" appears. Click "OK" to continue ■ If the following error messages are displayed at the end of the backup routine, ignore them. The backup is then completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
■ Remove the CD.	☐
■ Click Home and close the local service UI.	☐



If the backup creation to CD fails, please create the master backup to hard disk and transfer it to your service PC. On your service PC, burn the master backup to CD.

MRAWP Step: Backing up Dynamic Data at the MRAWP (MRC)






If you have a satellite console with a shared database connected to the main console, start with the action item:

- (→ MRWP Step: Backing up Dynamic Data at the MRWP (MRSC) / Page 33)

This item has to be performed first. If this step has already been completed on the MRSC, continue with the MRC.

Collecting MRAWP Configuration Data from Local Service

ACTION	MRAWP
■ Select Options --> Service --> Local Service.	☐
■ Enter the valid service key and click OK.	☐
■ Select Reports in the service platform.	☐

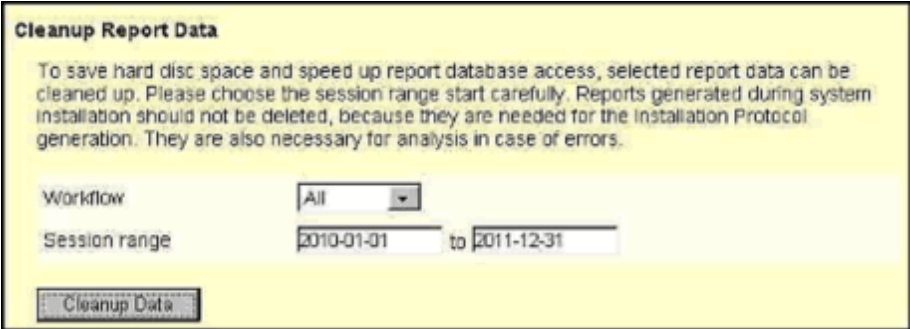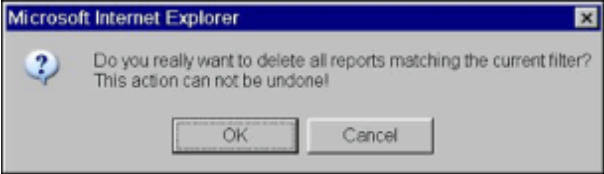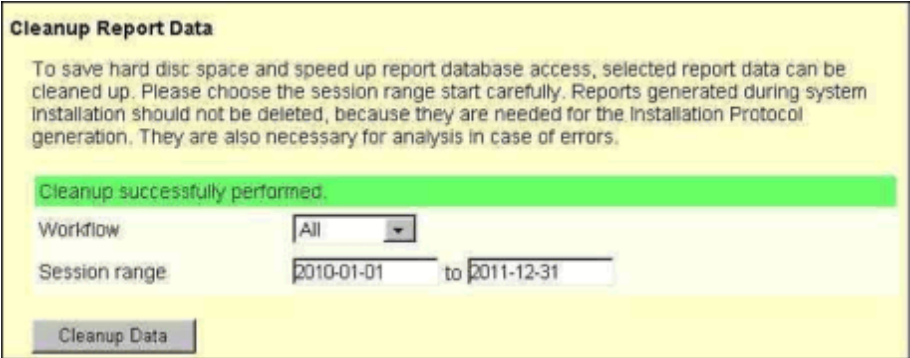ACTION	MRAWP
■ Purge service-related files (test tool files and tune-up, QA files older than one year). If the files are not purged, the size of the master backup created will be huge. - Select Reports -> Clean up - Under Cleanup Report Data select Workflow All - Session range 2010-01-01 to (today - 1 year) e.g. 2011-07-01 - You can also click in the fields and a calendar popup appears. <i>Fig. 33: Data Cleanup</i>  (1) Start date (2) End date - Press the Cleanup Data button. <i>Fig. 34: Confirmation popup</i>  - Confirm the popup with OK. - The following shows a successful file purge. <i>Fig. 35: Successful cleanup message</i> 	<input data-bbox="1353 275 1390 313" type="checkbox"/>
■ On the Report Files page, select Site Info Imprint	<input data-bbox="1353 1836 1390 1874" type="checkbox"/>

ACTION	MRAWP
Print out the Imprint report (generate a PDF file or write down the system options) Launch the update of the Imprint file (bottom right corner) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Imprint.pdf	☐
On the Report Files page, select Site Info Configuration	☐
- Print out the Configuration report (generate a PDF file or write down the system options) - DICOM networking (write down all DICOM network nodes and settings) - DICOM print devices (write down all DICOM print devices and settings) - DICOM HIS/RIS (write down all DICOM HIS/RIS nodes and settings) In order to generate a PDF file of the Configuration report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\Report.pdf	☐
On the Report Files page, select Site Info MR Configuration	☐
Print out the MR Configuration Report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MrReport.pdf.	☐
On the Report Files page, select Site Info CAN24 Configuration	☐
Print out the CAN24 Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\CAN24Report.pdf	☐
On the Report Files page, select Site Info MBox Configuration	☐
Print out the MBox Configuration report (generate a PDF file or write down the system options) In order to generate a PDF file of the Imprint report, press "show as PDF". Save the PDF file under C:\MedCom\MriSiteData\MBoxReport.pdf	☐
Click EXIT to close the Report Files popup window	☐



It is recommended to copy the PDF report files to an external medium to have the information for the configuration check in the new software version.

Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

Master Backup at the MRAWP (MRC)

ACTION	MRAWP
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Master Backup VD13 or VE11A or VE11B MRAWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK.	☐
Note: - Check the output log for errors. If any other error occurs apart from those mentioned below, the master backup has to be created again. - In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed" appears. Click "OK" to continue - If the following error messages are displayed at the end of the backup routine, ignore them. The backup is then completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed .. and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
Remove the CD.	☐
Click Home and close the local service UI.	☐



If the backup creation to CD fails, please create the master backup to hard disk and transfer it to your service PC. On your service PC, burn the master backup to CD.

Sending the Backup CDs to CV ME

Send the backup CDs to the following address:

Siemens Healthcare GmbH

HC DI CV ME OMM CPGMP

Attn.: Mr. G. Pfaendner

Allee am Roethelheimpark 2

D - 91052 Erlangen

"EVOLVE/CSSI Customer Data"

2.2.1.3 Perform remote backup for the CSSI update



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the requirement for full access. Also ask the customer to keep the system running in standby mode for the backup transfer.

Prerequisites

- System is connected to the SRS infrastructure via a VPN connection
- System is in **full remote access mode**
- The estimated time is about 30-40 min. for the remote backup, depending on the amount of backup data.

SRS Configuration Check

- Remote login to the system using the cRSP (common Remote Service Platform)
 - Check the RDIAG function on the **MRAWP** and, if applicable, on the **MRWP** as well
 - Test the FTP connection to the RDIAG server as follows:
 - Local Service -> File&Image Tools
 - Source: File
 - Action: Transfer to System
 - File Type: Any file
 - From source: SRS_ftp
- If a list from the directory /outgoing/mr (directory on the RDIAG server) is displayed, the connection and FTP access are OK



If the FTP test fails, check the SRS configuration; refer to (→ Final Check - to be performed for MRC and MRSC (if applicable) / Page 173)



If you do not have a satellite console with a shared database connected to the main console, the step (→ Auto Report setup and remote backup on MRWP / MR00-000.814.55) can be skipped.

Continue then with (→ MRAWP Step: Auto Report Setup and Remote Backup / Page 46).

MRWP Step: Auto Report Setup and Remote Backup



To create the backup, full access to the system is required. Please inform the customer and request full access. The time required for full access mode is about 30 minutes.

Setup of the CSSI Auto Report at the MRWP



The setup of the CSSI autoreport is only necessary if this function is not available on the system. Please check if the CSSI autoreports are enabled.

ACTION	MRWP
■ Login to the system via cRSP (common Remote Service Platform).	☐
■ Enter the valid service key and click OK .	☐
■ Select Auto Report in the service platform.	☐
■ Under Report , select the entry <code>TEMPLATE</code> .	☐
■ Create the new Auto Report with the name <code>CSSI_Trان_N4</code> Use the following values: - Active: Check this box - Name: <code>CSSI_Trان_N4</code> - Target: <code>SRS_mail</code> - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: <code>C:\MedCom\service\transfer*Numaris4*.ar</code> - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: <code>SRS_ftp</code> - Send when: At once	☐

ACTION

MRWP

Fig. 36: CSSI_Tran_N4


- Click **Save**.

☐

- Confirm the popup with **OK**.

☐

Fig. 37: Save OK

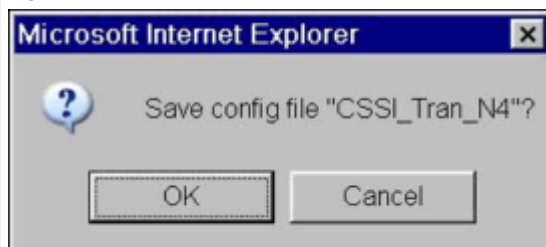

- Continue by clicking **OK** on the popup.

☐

Fig. 38: Success pop-up





- Under **Report**, select the entry **TEMPLATE**.




☐

ACTION	MRWP
■ Create the new Auto Report with the name CSSI_Trans_SecSe Use the following values: - Active: Check this box - Name: CSSI_Trans_SecSe - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*Security-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

Fig. 39: Autoreport CSSI_Trans_SecSe

■ Click Save.	☐
--	--------------------------

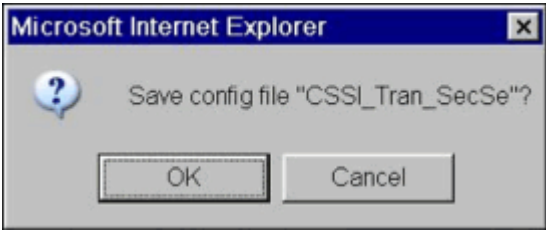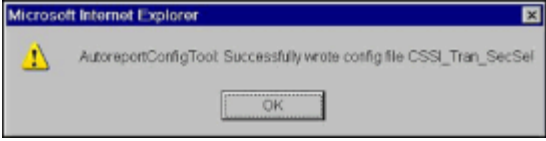ACTION	MRWP
Confirm the popup with OK. Fig. 40: Save OK 	☐
Continue by clicking OK on the popup. Fig. 41: Success pop-up 	☐
Under Report, select the entry TEMPLATE.	☐
- Create the new Auto Report with the name <code>CSSI_Tran_SW02</code> Use the following values: - Active: Check this box - Name: <code>CSSI_Tran_SW02</code> - Target: <code>SRS_mail</code> - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: <code>C:\MedCom\service\transfer*SW-Settings*.ar</code> - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: <code>SRS_ftp</code> - Send when: At once	☐

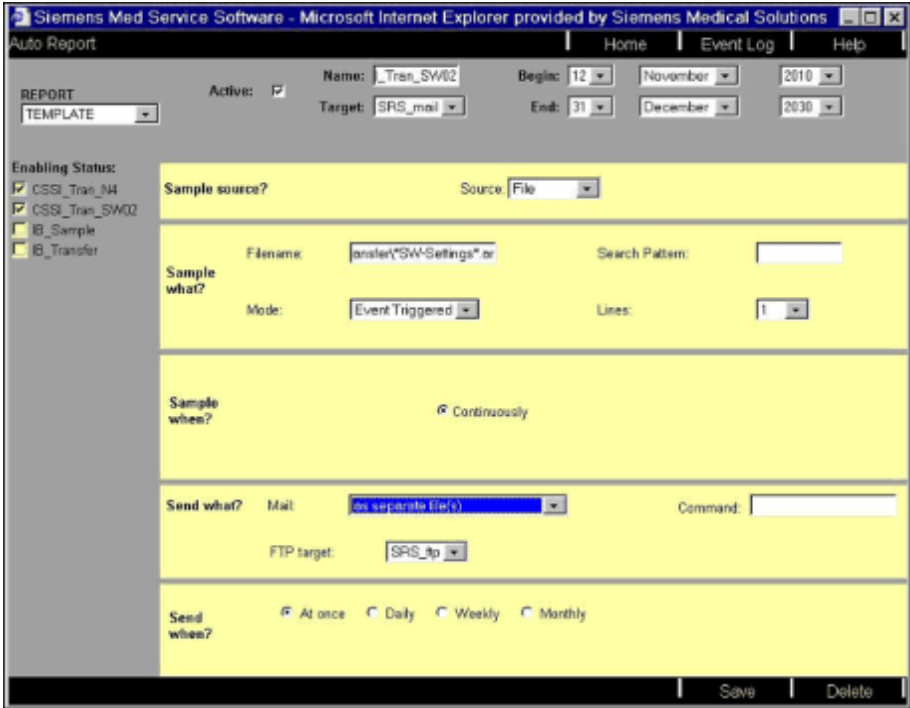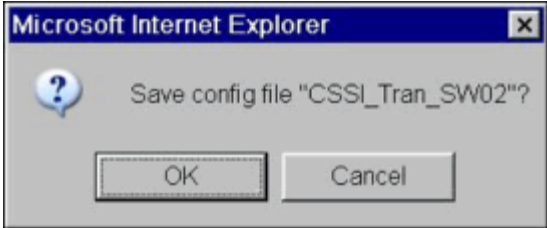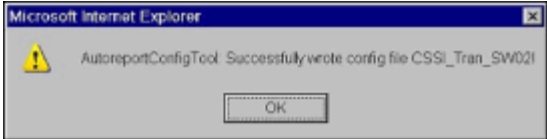ACTION		MRWP
Fig. 42: Autoreport CSSI_Tran_SW02 		
■ Click Save .		☐
■ Confirm the popup with OK .		☐
Fig. 43: Save OK 		
■ Continue by clicking OK on the popup.		☐
Fig. 44: Success pop-up 		
■ Check the newly created Autoreports to see if they are configured correctly and are activated.		☐
■ Click Home to leave the Auto Report section.		☐

Performing the Backup at the MRWP



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the requirement for full access. Also ask the customer to keep the system running in standby mode for the backup transfer.

ACTION	MRWP
Select Backup&Restore in the service platform.	☐
- Perform the backup in local service. The backup has to be performed for the following packages: Numaris4, SW settings02, and Security Settings. Select Backup in the Command menu, select [-E-] FIXED in the Drives menu. - Select SW Settings02 in the Package menu, start the backup by clicking OK. - Select Security Settings in the Package menu, start the backup by clicking OK. - Select Numaris4 in the Package menu, start the backup by clicking OK. Check for errors using the log window.	☐
Click Home to leave Backup&Restore	☐

Move Backups into the transfer directory MRWP

ACTION	MRWP
Select Utilities in the service platform.	☐
Under Source, select the entry Escape to OS in the drop-down menu	☐
Under Command, keep the default entry NT Command Interpreter.	☐
Under Parameters, type in the following: <code>move /y E:*.ar C:\MedCom\service\transfer\</code>	☐

ACTION	MRWP	
Click Go. The three backups that have been created are now moved to the directory C:\Med-Com\service\transfer. Fig. 45: Move output  <tr><td> Click Home to leave Utilities Check the Auto Reports to see whether they are active or not (see (→ Checking active CSSI Auto Reports / MR00-000.814.55)) Close the service UI in order to close the connection. </td><td> ☐ ☐ ☐ </td></tr>	Click Home to leave Utilities Check the Auto Reports to see whether they are active or not (see (→ Checking active CSSI Auto Reports / MR00-000.814.55)) Close the service UI in order to close the connection.	☐ ☐ ☐
Click Home to leave Utilities Check the Auto Reports to see whether they are active or not (see (→ Checking active CSSI Auto Reports / MR00-000.814.55)) Close the service UI in order to close the connection.	☐ ☐ ☐	

MRAWP Step: Auto Report Setup and Remote Backup

Setup of the CSSI Auto Report at the MRAWP





The setup of the CSSI autoreport is only necessary if this function is not available on the system. Please check if the CSSI autoreports are enabled.

ACTION	MRAWP
Login to the system via cRSP (common Remote Service Platform).	☐
Select Auto Report in the service platform.	☐
Under Report , select the entry <code>TEMPLATE</code> .	☐

ACTION	MRAWP
■ Create the new Auto Report with the name CSSI_Trان_N4 Use the following values: - Active: Check this box - Name: CSSI_Trان_N4 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*Numaris4*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐

Fig. 46: CSSI_Trان_N4

■ Click Save.	☐
--	--------------------------

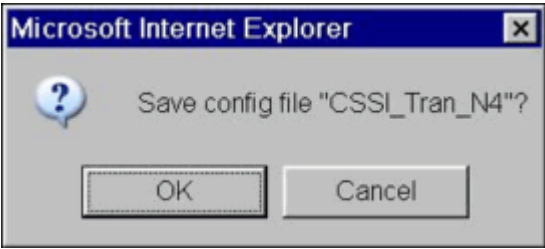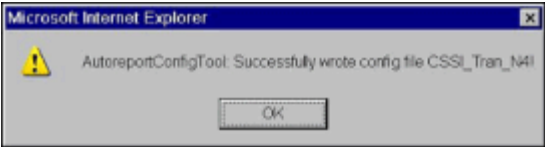ACTION	MRAWP
Confirm the popup with OK. Fig. 47: Save OK 	☐
Continue by clicking OK on the popup. Fig. 48: Success pop-up 	☐
Under Report, select the entry <code>TEMPLATE</code>.	☐
- Create the new Auto Report with the name <code>CSSI_Trans_SecSe</code> Use the following values: - Active: Check this box - Name: <code>CSSI_Trans_SecSe</code> - Target: <code>SRS_mail</code> - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: <code>C:\MedCom\service\transfer*Security-Settings*.ar</code> - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: <code>SRS_ftp</code> - Send when: At once	☐

ACTION

MRAWP

Fig. 49: Autoreport CSSI_Tran_SecSe


- Click **Save**.

☐

- Confirm the popup with **OK**.

☐


Fig. 50: Save OK

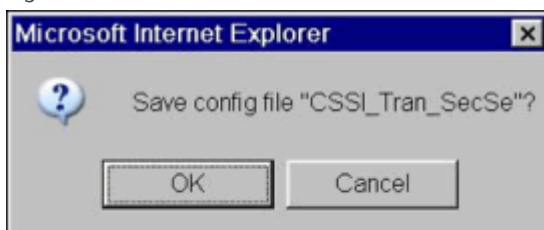

- Continue by clicking **OK** on the popup.


☐



Fig. 51: Success pop-up

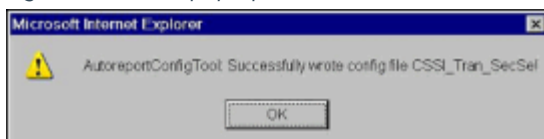

- Under **Report**, select the entry **TEMPLATE**.

☐

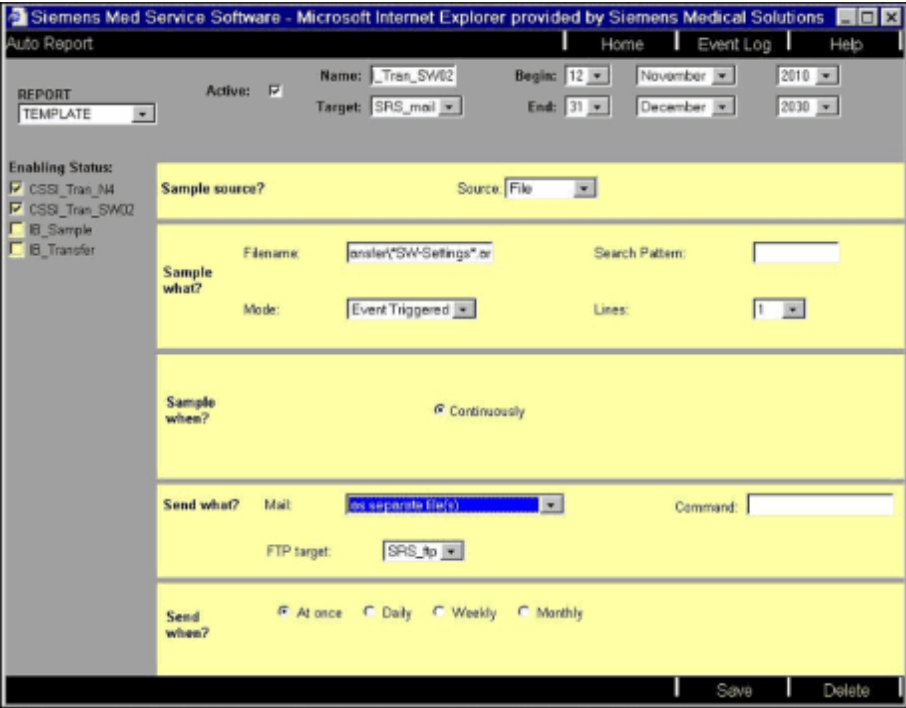ACTION	MRAWP
■ Create the new Auto Report with the name CSSI_Trans_SW02 Use the following values: - Active: Check this box - Name: CSSI_Trans_SW02 - Target: SRS_mail - Begin: No change required (date is today) - End: 31 December 2030 - Source: File - File name: C:\MedCom\service\transfer*SW-Settings*.ar - Mode: Event triggered - Sample when: Continuously - Mail: As separate file(s) - FTP target: SRS_ftp - Send when: At once	☐
Fig. 52: Autoreport CSSI_Trans_SW02 	☐
■ Click Save.	☐


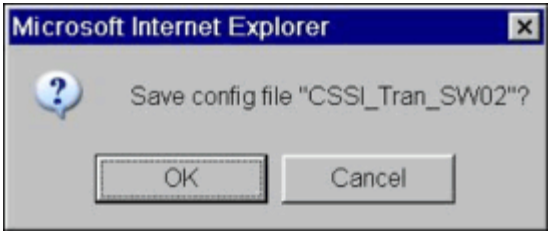
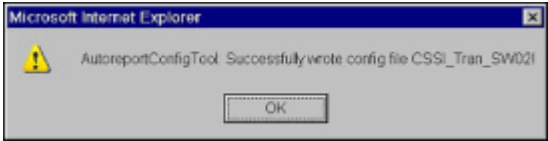
ACTION	MRAWP
Confirm the popup with OK. <i>Fig. 53: Save OK</i> 	☐
Continue by clicking OK on the popup. <i>Fig. 54: Success pop-up</i> 	☐
Check the newly created Autoreports to see if they are configured correctly and are activated.	☐
Click Home to leave the Auto Report section.	☐

Performing the Backup at the MRAWP



Before the remote backup can be performed, you must inform the customer about the duration of the backup and the requirement for full access. Also ask the customer to keep the system running in standby mode for the backup transfer.

ACTION	MRAWP
Select Reports in the service platform.	☐

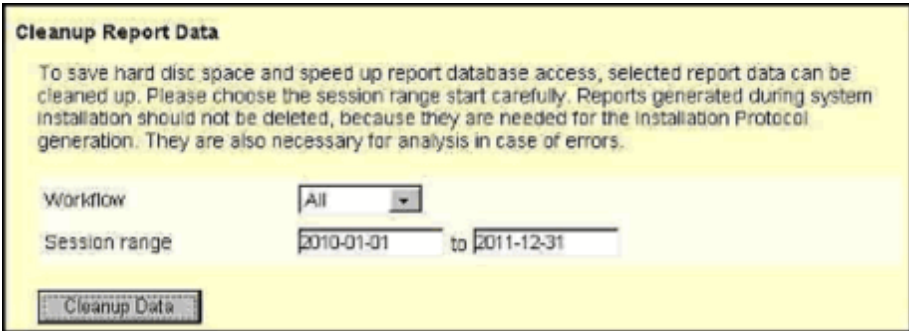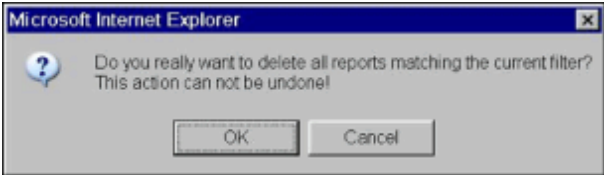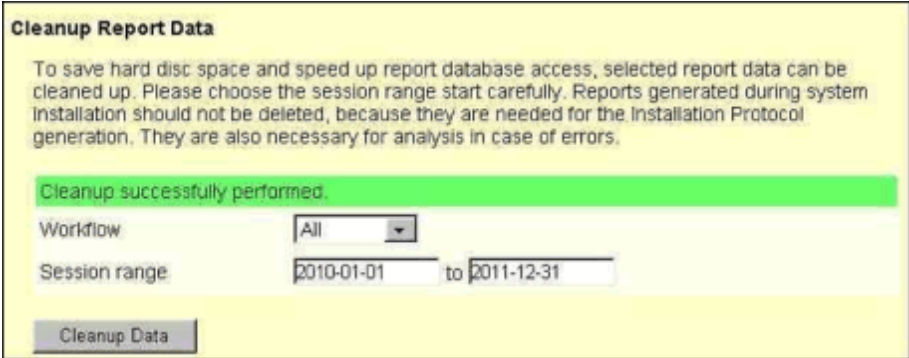ACTION	MRAWP
■ Purge service-related files (test tool files and tune-up, QA files older than one year). If the files are not purged, the size of the master backup created will be huge. - Select Reports -> Clean up - Under Cleanup Report Data select Workflow All Session range 2010-01-01 to (today - 1 year) e.g. 2011-07-01 You can also click in the fields and a calendar popup appears. Fig. 55: Data Cleanup  (1) Start date (2) End date - Press the Cleanup Data button. Fig. 56: Confirmation popup  Confirm the popup with OK. - The following shows a successful file purge. Fig. 57: Successful cleanup message 	☐
■ Click EXIT to close the Report Files popup window	☐
■ Select Backup&Restore in the service platform.	☐

ACTION	MRAWP
- Perform the backup in local service. The backup has to be performed for the following packages: Numaris4, SW settings02, and Security Settings. Select Backup in the Command menu, select [-E-] FIXED in the Drives menu. - Select SW Settings02 in the Package menu, start the backup by clicking OK. - Select Security Settings in the Package menu, start the backup by clicking OK. - Select Numaris4 in the Package menu, start the backup by clicking OK. The duration of the Numaris4 backup is about 15-20 min. depending on the amount of data. Check for errors using the log window.	☐
Click Home to leave Backup&Restore	☐

Move Backups into the transfer directory

ACTION	MRC
Select Utilities in the service platform.	☐
Under Source, select the entry Escape to OS in the drop-down menu	☐
Under Command, keep the default entry NT Command Interpreter.	☐
Under Parameters, type in the following: <code>move /y E:*.ar C:\MedCom\service\transfer\</code>	☐

ACTION	MRC
■ Click Go. The three backups that have been created are now moved to the directory C:\Med-Com\service\transfer. Fig. 58: Move output	☐
■ Click Home to leave Utilities	☐

Checking active CSSI Auto Reports

ACTION	MRAWP/MRWP
■ Select File & Image Tools in the service platform.	☐
■ Under Source, select the entry Transfer List from the drop-down menu.	☐
■ Under Action, select the entry Show/Delete in the drop-down menu.	☐
■ Now the Transfer List is shown. Check to see if the FTP transfer for the backup files has been started. To do this, select the FTP entries in the transfer manager list.	☐

ACTION

MRAWP/MRWP

■ For the Numaris4 backup file.

Fig. 59: Transfer List active autoreport Numaris4

The screenshot shows the 'File & Image Tools' window with the 'Transfer List' active autoreport for Numaris4. The window displays a list of transfer entries and detailed transfer information for a specific entry.

Entry time	Source file name(s)	Transfer type
Mon Nov 15 08:12:40 2010	C:/MedCom/service/Backup/output/	Mail
Mon Nov 15 08:12:37 2010	C:/MedCom/service/transfer/CSF/	FTP system-target
Mon Nov 15 08:12:36 2010	C:/MedCom/service/Backup/output/	Mail
Mon Nov 15 08:12:35 2010	C:/MedCom/service/transfer/CSF/	FTP system-target
Mon Nov 15 08:12:33 2010	C:/MedCom/service/Backup/output/	Mail

Transfer Manager queue ID: 6
Entry time: Mon Nov 15 08:12:33 2010
Source file name(s): C:/MedCom/service/transfer/CSF/PSK/YRA1_40478_Numaris4_20101115_074612_ek.er
Target: sap00670
Transfer type: FTP system-target
Transfer status: Transferring
Transfer status details: transfer thread started
Abort flag: Queued
Total bytes to transfer: 59532160
Actually transferred bytes: 1116160
Retry limit: 5
Retry count: 1
Transfer time: immediately

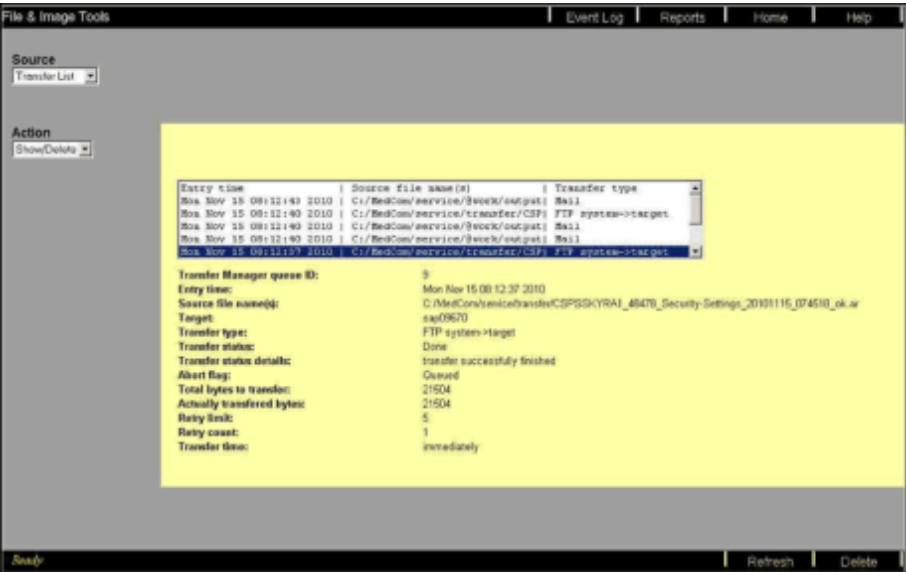
■ For the Software02 backup file.

Fig. 60: Transfer List active autoreport Software02

The screenshot shows the 'File & Image Tools' window with the 'Transfer List' active autoreport for Software02. The window displays a list of transfer entries and detailed transfer information for a specific entry.

Entry time	Source file name(s)	Transfer type
Mon Nov 15 08:12:43 2010	C:/MedCom/service/Backup/output/	Mail
Mon Nov 15 08:12:40 2010	C:/MedCom/service/transfer/CSF/	FTP system-target
Mon Nov 15 08:12:40 2010	C:/MedCom/service/Backup/output/	Mail
Mon Nov 15 08:12:40 2010	C:/MedCom/service/transfer/CSF/	FTP system-target
Mon Nov 15 08:12:37 2010	C:/MedCom/service/Backup/output/	Mail

Transfer Manager queue ID: 12
Entry time: Mon Nov 15 08:12:40 2010
Source file name(s): C:/MedCom/service/transfer/CSF/PSK/YRA1_40478_SW-Settings02_20101115_074422_ek.er
Target: sap00670
Transfer type: FTP system-target
Transfer status: Transferring
Transfer status details: transfer thread started
Abort flag: Queued
Total bytes to transfer: 2999000
Actually transferred bytes: 1085440
Retry limit: 5
Retry count: 1
Transfer time: immediately

ACTION	MRAWP/MRWP
■ For the Security Settings backup file. Fig. 61: Transfer List active autoreport SecuritySettings 	☐
■ Click Home to leave Utilities	☐
■ Close the local service UI.	☐



Please ask the customer not to switch off the overall system. The system should be switched to standby.

3.1 CSSI Upgrade Second Steps for VE11A/VE11B




After receiving the CSSI kit for the software upgrade to VE11A/VE11B, the following steps have to be performed.

3.1.1 Work Steps

3.1.1.1 Overview

Perform the CSSI upgrade

Fig. 62: CSSI Final Steps Flowchart

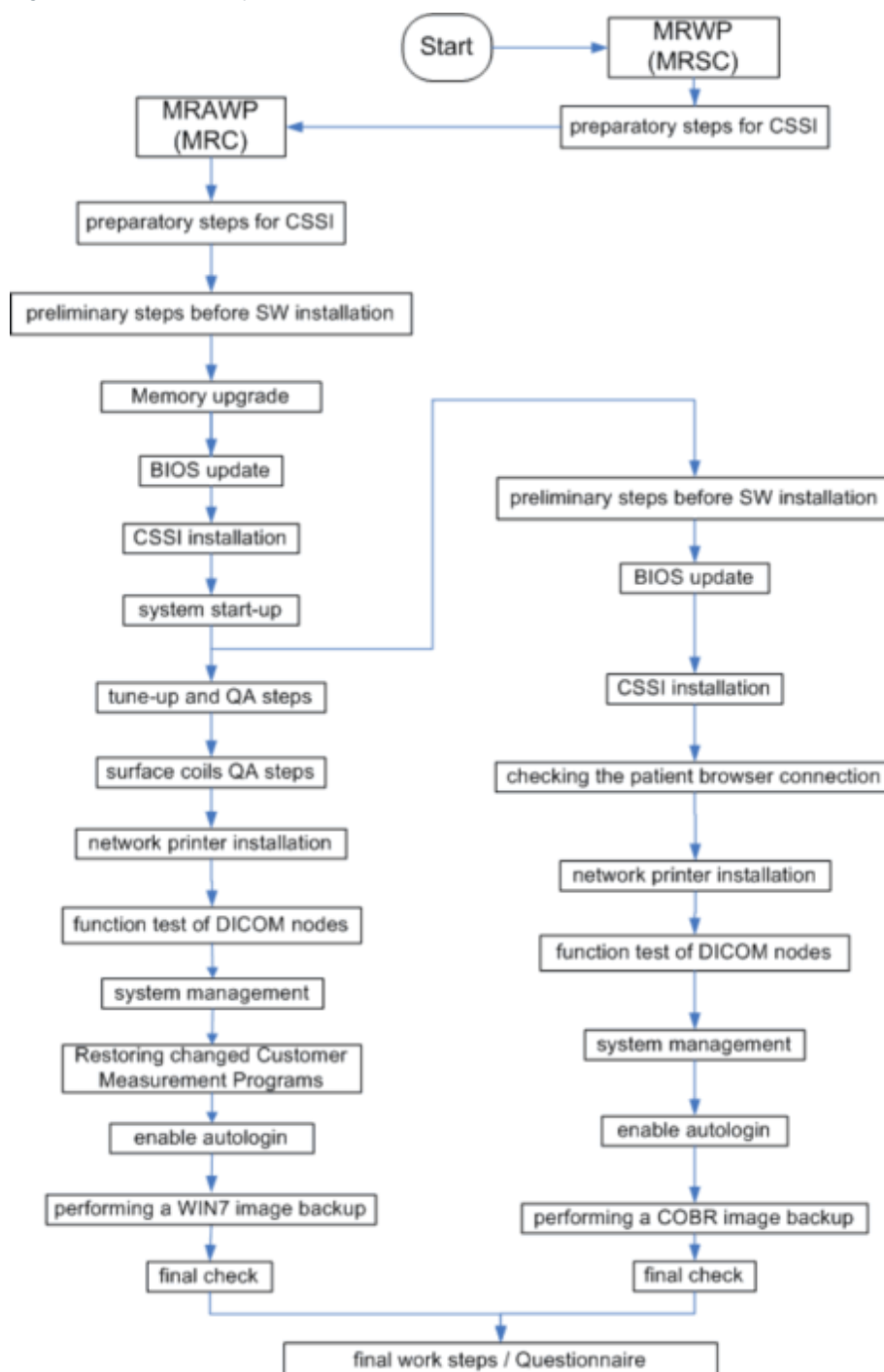

3.1.1.2 MRWP (MRSC) preparatory Steps for the CSSI Upgrade with the old Software Version



If you have a satellite console with a shared database connected to the main console, perform the following steps first.

If no MRWP is present, skip the following steps and start at (→ MRAWP (MRC) preparatory Steps for CSSI Update with the old Software Version / Page 60)



Please check if any local network printer is configured on MRWP and clarify with the customer representative the availability of installation material (Software, Installation Manual,..) if necessary. Please note the configuration data.

Performing a Master Backup



Please check if any local or network printer is configured on the MRC and clarify with the customer representative the availability of installation material (Software, Installation Manual,..) if necessary. Please note the configuration data.



If required (configuration is changed) the step (→ Collecting MRWP Configuration Data from Local Service / Page 10) can be done too.

ACTION	MRWP
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Master Backup VD13A or VD13B or VD13C or VD13D or VD13E MRWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK.	☐
Note: - In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed," appears. Click "OK" to continue - If the following error messages are displayed at the end of the backup routine, ignore them. The backup is completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
Remove the CD.	☐

MNP Deregistration on the MRWP



In order to avoid problems with the MNP registration with the new Numaris version it is recommended to perform the deregistration with the old Numaris version.

ACTION	MRWP
■ In the local Service UI, select Configuration .	☐
■ Under System Mgmt. Agent Controls , set the following configuration: Remove the mark from the checkboxes "Go Active" and "Go Online".	☐
■ Click Go .	☐
The deregistration of the MNP is done. This will take about 5 minutes.	
■ Click Finish .	☐
■ Click Home and close the local service UI.	☐
■ Shut down the MRWP.	☐
■ Switch off the MRWP.	☐

3.1.1.3 MRAWP (MRC) preparatory Steps for CSSI Update with the old Software Version



If you have a satellite console with a shared database connected to the main console, start with the action item:

(→ MRWP (MRSC) preparatory Steps for the CSSI Upgrade with the old Software Version / Page 58)

This has to be done first. If this step has already been completed on the MRWP, continue with the MRAWP.



Please check if any local network printer is configured on MRAWP and clarify with the customer representative the availability of installation material (Software, Installation Manual,..) if necessary. Please note the configuration data.


Performing a Master Backup on the MRC, for changed Customer Protocols



Please check if any local or network printer is configured on the MRC and clarify with the customer representative the availability of installation material (Software, Installation Manual,..) if necessary. Please note the configuration data.



If required (configuration is changed) the step (→ Collecting MRAWP Configuration Data from Local Service / Page 12) can be done too.

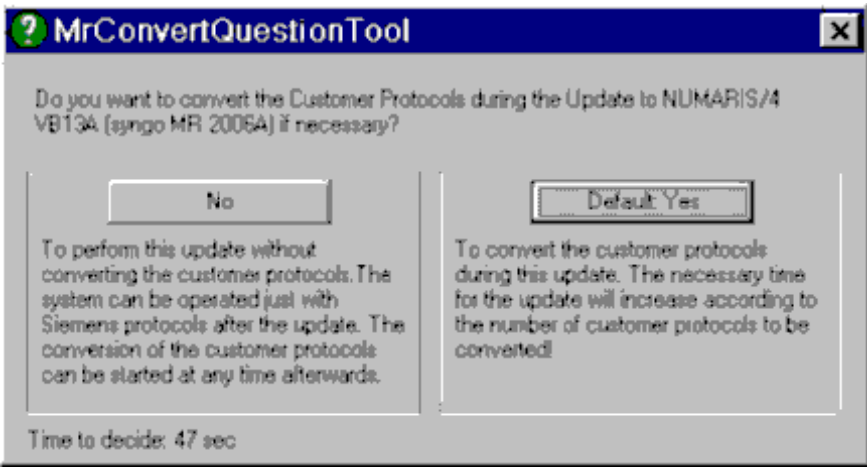ACTION	MRAWP
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Master Backup VD13A or VD13B or VD13C or VD13D or VD13E MRAWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK. During the Master Backup, the following popup appears: Fig. 63:  Press the "Yes" button. At this point, this popup will not start the protocol conversion if "yes" was selected. This popup determines the default focus ("Yes" or "NO") of the MrConvertQuestion-Tool when restoring the Master Backup.	<input data-bbox="1380 275 1420 309" type="checkbox"/>
Note: - In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed," appears. Click "OK" to continue - If the following error messages are displayed at the end of the backup routine, ignore them. The backup is completed correctly. - Error Execute starting components... 3 rti.exe runComp an_serv- er rti->runComp(an_server, Pid, Container) failed and/or - Error .. Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed .. and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ..	
Remove the CD.	<input data-bbox="1380 1917 1420 1951" type="checkbox"/>

MNP Deregistration on the MRAWP



In order to avoid problems with the MNP registration with the new Numaris version it is recommended to perform the deregistration with the old Numaris version.

ACTION	MRAWP
In the local Service UI, select Configuration.	☐
Under System Mgmt. Agent Controls, set the following configuration: Remove the mark from the checkboxes "Go Active" and "Go Online".	☐
Click Go.	☐
The deregistration of the MNP is done. This will take about 5 minutes.	
Click Finish.	☐
Click Home.	☐

3.1.1.4 MRAWP (MRC) Preliminary Steps before installing the Software

ACTION	MRAWP
Be sure that the network cables at the back of the console are connected. - MRC: Two network cables at the back of the computer.	☐
Be sure that only the monitor and the PDD are connected. Disconnect any additionally connected USB customer device (like USB harddisks). Keep only the USB mouse, the USB keyboard and powerloading stations for ECG connected.	☐
If applicable, check whether the local connected printer (connected to the parallel port) is switched ON. During the installation, the printer is automatically detected and the correct printer driver is installed.	☐
Disable the Autologin function in the local service UI.	☐
Perform a shutdown and switch of the system on the alarmbox.	☐
Remove the power line connector from the MARS in the Electronic Power Cabinet. (The reason for it is that the MARS installation can run into a problem if it is not done)	☐

3.1.1.5 MRAWP (MRC) Memory upgrade to 32 GByte




Use the ESD kit when performing the upgrade.

Only valid for MRAWP HPz420 host with

- partnumber **10864002** (see Siemens partnumber sticker)
- or
- only valid for Brazil: partnumber **10928155** (see Siemens partnumber sticker)

Fig. 64: Example of partnumber sticker (Model No.)

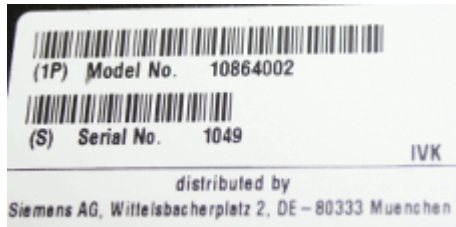

Disassemble HPz420

- Remove all connected cables
- Use ESD outfit.
- Open the casing
 - Remove the two lock screws on the back of the casing.


Fig. 65: HP Z420 rear view



- (1) Lock screw
- (2) Lock screw
- Lay the device on its side on a stable, flat, and clean surface.


- Pull up on the handle (1), slide the cover 1/2-inch toward the rear of the system (2), then rotate the cover off the chassis (3).

Fig. 66: HP Z420 side panel

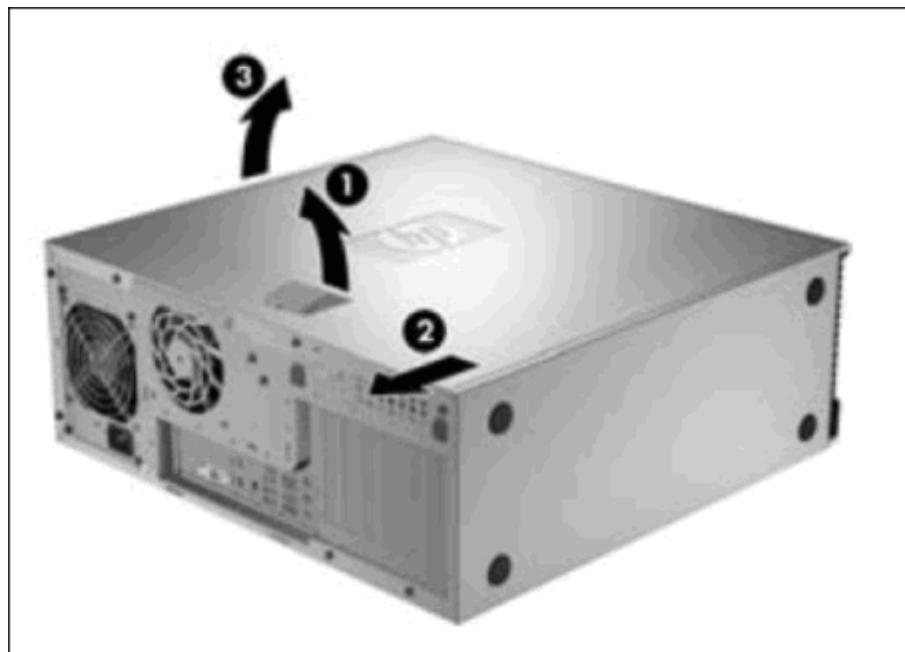

- Removing the fan from the second memory location


Fig. 67: HPz420 inside

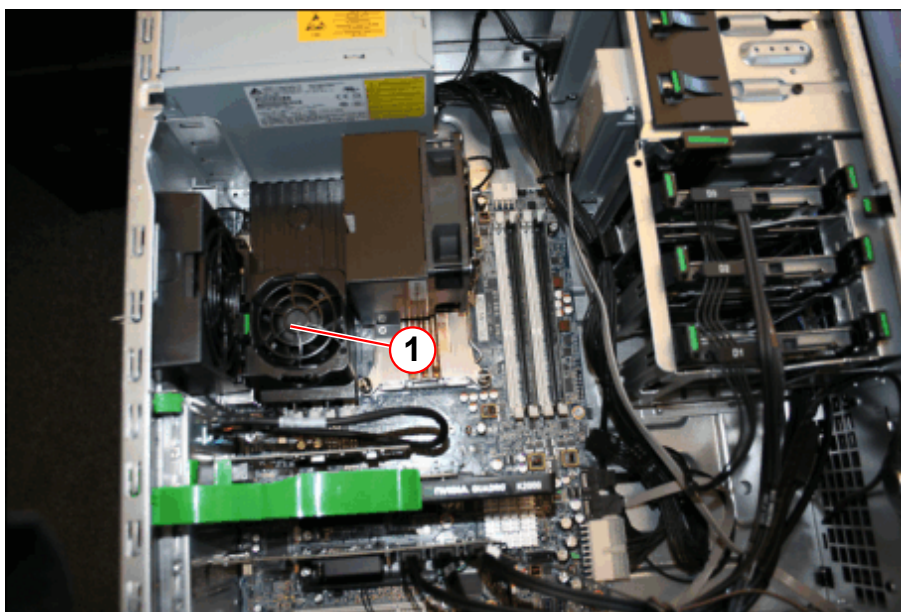

(1) Fan to remove


- Press the fan release latch as shown in the following figure, and guide fan out of the mounting location..

Fig. 68: Fan remove

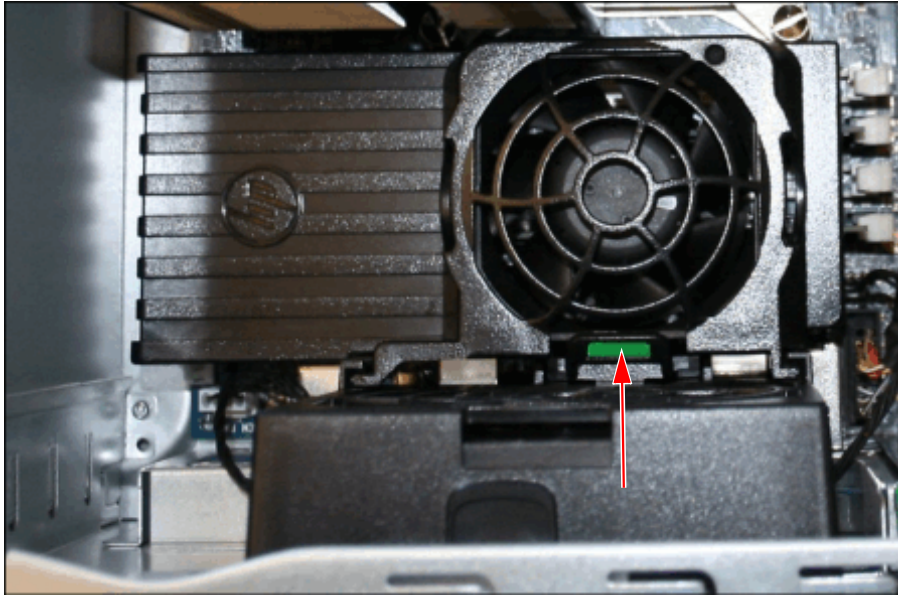

- Place the removed fan on the power supply (see following figure).

Fig. 69: HPz420 with removed Fan

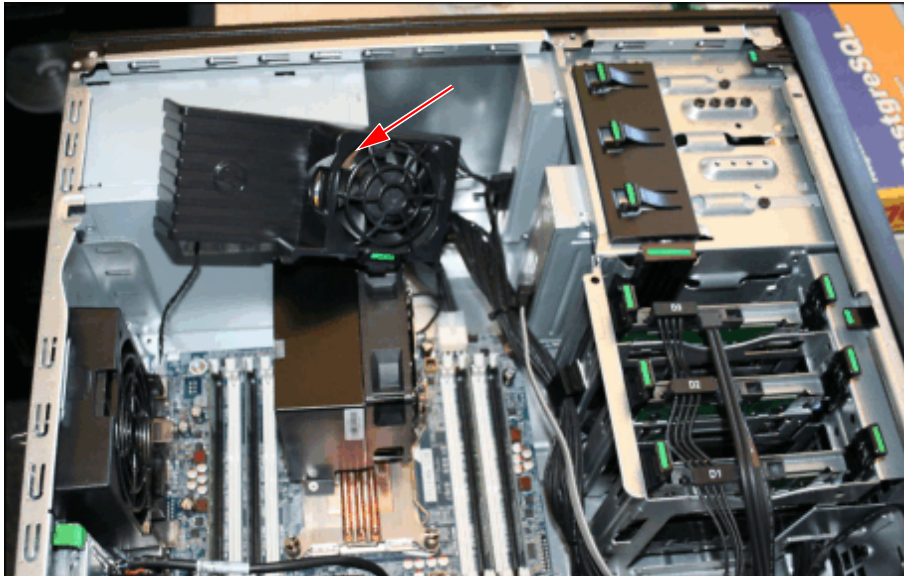

RAM installation for HPz420

- Install the additional 4 GB memory modules in the empty module locations (white sockets).




Only the following memory module can be used:

10928151 MEMORYMOD 4GB DDR3-1600 UNB (Siemens part number and description)

HP description: HP 4GB (1x4GB) DDR3-1600 ECC or DDR3-1866 UNB (both are usable)

Fig. 70: HPz420 memory module locations

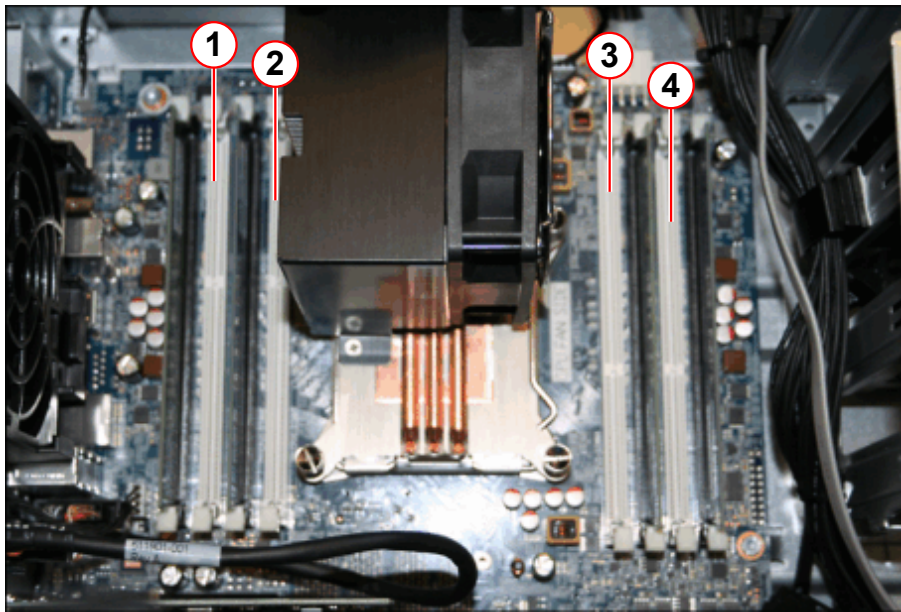

- (1) Place 4 GB memory module white empty DIMM socket
- (2) Place 4 GB memory module white empty DIMM socket
- (3) Place 4 GB memory module white empty DIMM socket
- (4) Place 4 GB memory module white empty DIMM socket


Insert the 4 GB memory modules in the white empty DIMM socket.

When installing memory modules, make sure the modules align properly with the memory socket. There should be keys (small indents) on your memory modules that fit the keys in the memory socket. DDR modules and sockets have only one key which is near the center of the module/socket.

Procedure to install a DIMM memory module

- Push gently outward on the DIMM socket levers as shown in the following illustration. (The example shows a HP Z600 Workstation.)

Fig. 71: Opening DIMM socket levers

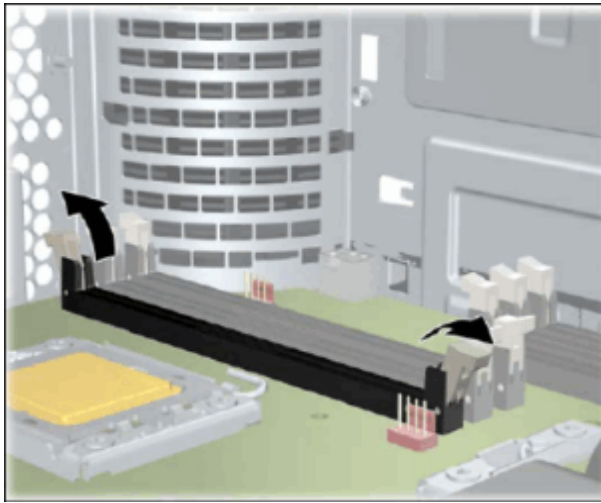


- Align the DIMM connector key with the DIMM socket key, and then seat the DIMM firmly in the socket (1) as shown in the following illustration.



DIMMs and their sockets are keyed for proper installation. To prevent socket or DIMM damage, align these guides properly when installing DIMMs.

Fig. 72: Seating the DIMM

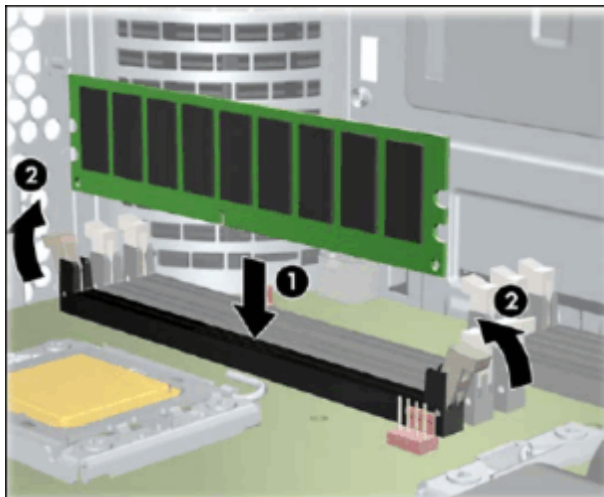

- Secure the socket levers (2).

Reassemble HPz420

- Installing the fan
 - Ensure that the latch on the Fan snaps in.
- Closing the casing
 - Lay the device on its side on a stable, flat, and clean surface.

- Hook the side cover into the cut-outs in the upper and lower guide rails. Pay attention to ensure that the lugs of the lock fit into the associated retainers on the casing.
- Slide in the side cover until you hear it engage in the clip.
- Secure the two lock screws of the side cover on the back of the casing.
- Relabel the host HPz420 using the new part number (10928295 or for Brazil only 10928294) sticker and new Rev. sticker..
- Return the device to the place where you require to use it.
- Reconnect any disconnected cables (power cord, cables to external devices, etc.).


ACTION	MRAWP
■ Memory installed	☐
■ New part number label placed on host.	☐

Memory check

- Check recognized memory size
 - Be sure that the power line connector is removed from the MARS.
 - Switch on the system.
 - Press F1 to skip the memory check from the BIOS.
 - Login as user administrator
 - On the desktop, select the icon "Computer"
 - Press the right mouse button and select `Properties` in the menu displayed.

- In the displayed property window, check the value for Installed memory (RAM) it should be:
 - 32.0 GB

Fig. 73: System properties

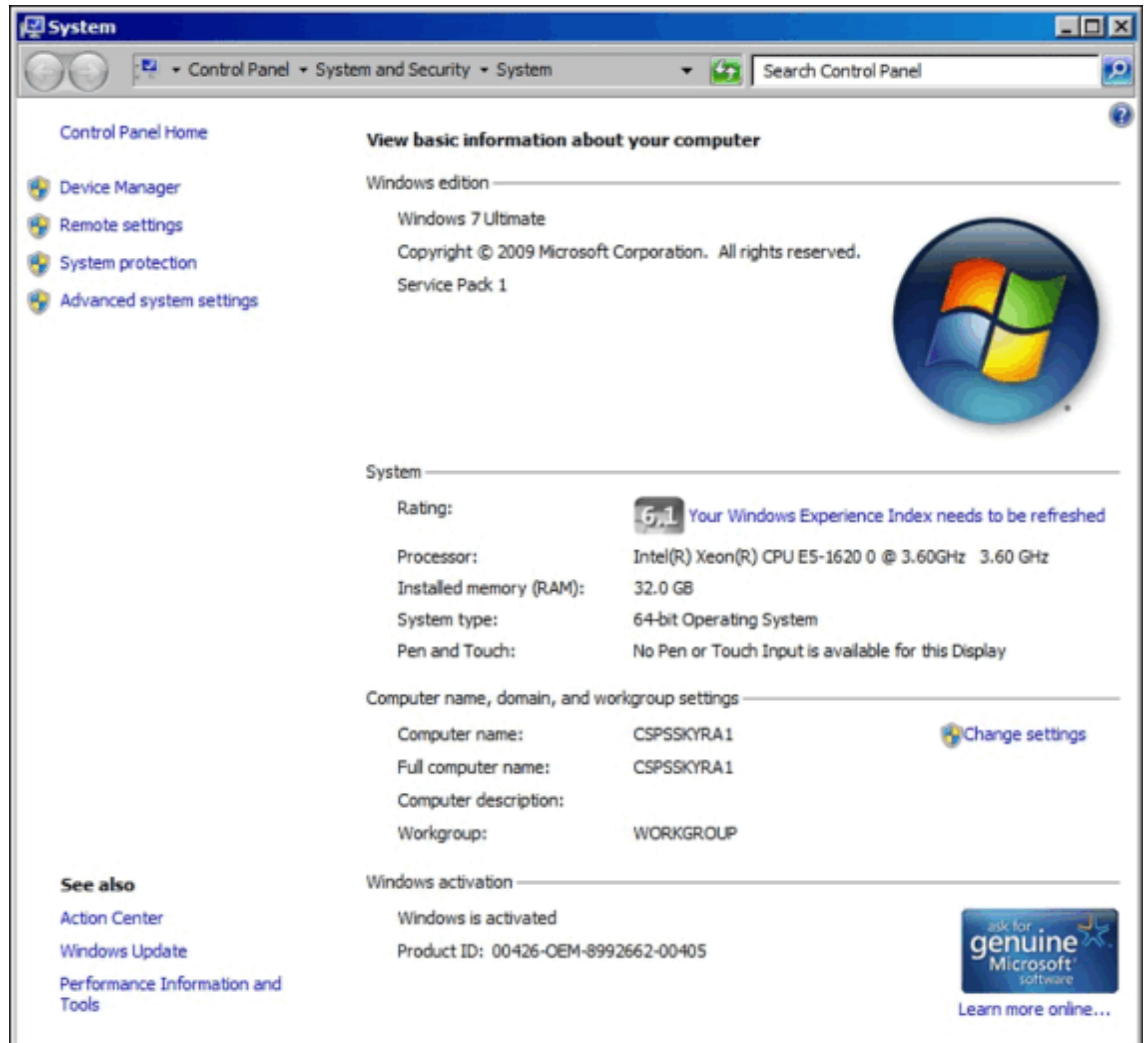


Is this not the case the memory modules are not correct installed.

- Close System properties popup window.

3.1.1.6 MRAWP (MRC) BIOS update



ACTION	MRAWP
■ Insert the first DVD of the delivered CSS image set labeled COBR Recovery Image MR system software VE11A or VE11B MR system serial no: <x> AWP into the CD/DVD burner drive (S:I).	☐

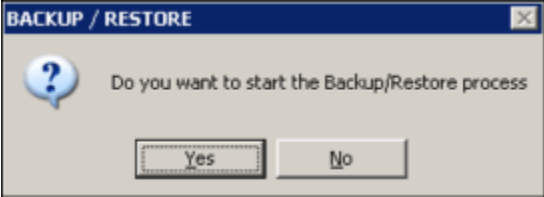ACTION	MRAWP
Insert the delivered CD labeled with "MR-VD BIOS CD VE11.." into the host CD/DVD-ROM Drive (R:\)	☐
Invoke a command shell (Start ⇒ All Programs ⇒ Accessories ⇒ Command Prompt)	☐
In the command shell, type in: S:\Start.bat	☐
Confirm the Installation popup with "Yes" <i>Fig. 74: Backup/Restore popup</i> 	☐
The running script performs a reboot at the end. This may take some time. Please wait!	☐
The BIOS CD will be booted automatically.	☐
The BIOS update runs automatically, no user action is necessary, and ends when the following message is displayed: <pre> BIOS Installation Done - System Ready for Software Installation Please remove BIOS CD ! Please Turn Off Power or Reset System ! C:\> </pre>	☐
Remove the delivered BIOS CD from the CD/DVD drive (R:\)	☐
Remove the first DVD of the delivered CSS image set labeled COBR Recovery Image MR system software VE11A or VE11B MR system serial no: <x> AWP from the CD/DVD burner drive (S:\).	☐

ACTION	MRAWP
Insert the first DVD of the delivered CSSI image set labeled COBR Recovery Image MR system software VE11A or VE11B MR system serial no: <x> AWP into the CD/DVD-ROM drive (R:).	☐
Press key combination <CTRL> + <ALT> + (on a German keyboard use <Strg> + <Alt> + <Entf>) to reboot and continue with "CSSI installation MRC".	☐

3.1.1.7 MRAWP (MRC) CSSI Installation

ACTION	MRAWP
A reboot is performed and the COBR tool boots.	☐
Note: If you get the following error message "Operating system not found", please remove the first DVD of the delivered CSS image set from the CD/DVD drive (R:) and insert it into the CD/DVD burner drive (S:) and perform a reboot.	
Wait until the COBR tool window is displayed.	☐
Fig. 75: COBR tool window 	

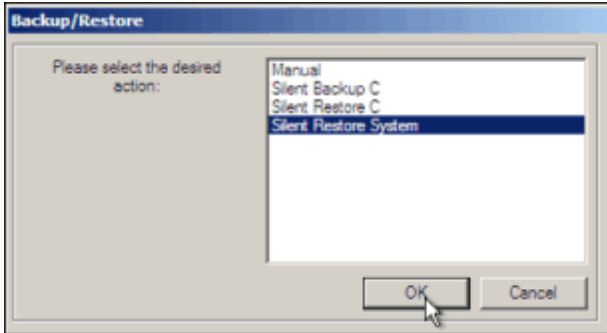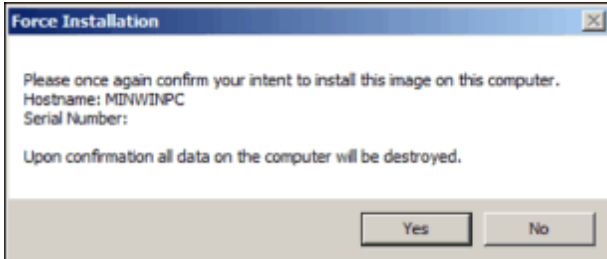
ACTION	MRAWP
In the Backup/Restore popup, select Silent Restore System and press OK Fig. 76: Backup/Restore popup Silent Restore System 	☐
Confirm the Force Installation popup (example is shown below) with Yes. Fig. 77: Force Installation popup  This Force Installation popup is displayed twice. Confirm it each time with Yes.	☐
The system restore process starts now. The following shows an example of the screen display. Fig. 78: COBR restore process start 	☐

ACTION

MRAWP

- First all image files from the DVD for the system disk (disk1) are copied from the DVD to the hard disk.

Fig. 79: Copy image files

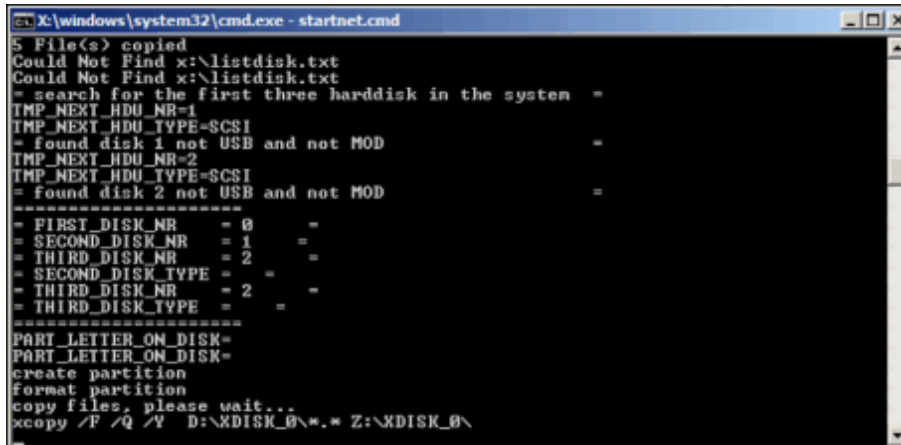

```

C:\windows\system32\cmd.exe - startnet.cmd
5 File(s) copied
Could Not Find x:\listdisk.txt
Could Not Find x:\listdisk.txt
- search for the first three harddisk in the system =
TMP_NEXT_HDU_NR=1
TMP_NEXT_HDU_TYPE=$CSI
- found disk 1 not USB and not MOD =
TMP_NEXT_HDU_NR=2
TMP_NEXT_HDU_TYPE=$CSI
- found disk 2 not USB and not MOD =
=====
- FIRST_DISK_NR = 0 =
- SECOND_DISK_NR = 1 =
- THIRD_DISK_NR = 2 =
- SECOND_DISK_TYPE = =
- THIRD_DISK_NR = 2 =
- THIRD_DISK_TYPE = =
=====
PART_LETTER_ON_DISK=
PART_LETTER_ON_DISK=
create partition
format partition
copy files, please wait...
xcopy /F /Q /Y D:\XDISK_0\*. * Z:\XDISK_0\


```

After the files are copied, the DVD is automatically ejected.

- Remove the ejected medium from the CD/DVD-ROM drive, insert the next DVD from the system image DVD set, and close the drive.

- Confirm the request popup with OK.


Fig. 80: Next DVD request popup

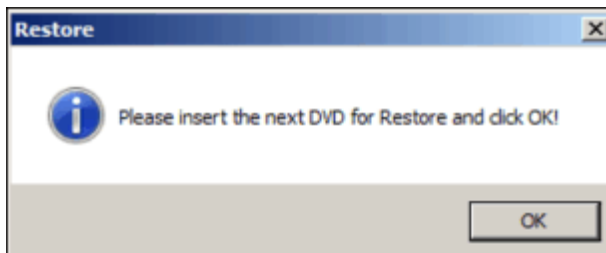

The copy procedure continues. The next DVD is requested automatically until all required data are copied. The copy and the confirmation of the popup (→ Fig. 80 Page 73) have to be repeated for all available CSSI DVDs.

- After the last CSSI DVD has been copied, the Image Restore starts



Fig. 81: Image restore progress output

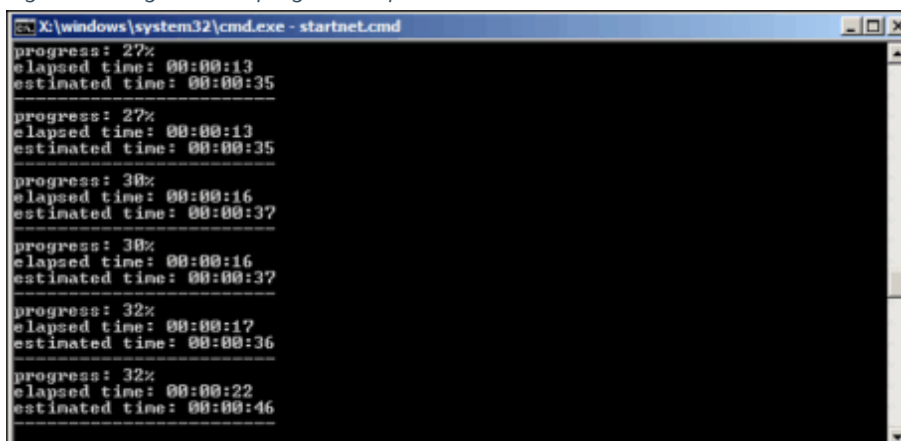

```

C:\windows\system32\cmd.exe - startnet.cmd
progress: 27%
elapsed time: 00:00:13
estimated time: 00:00:35
-----
progress: 27%
elapsed time: 00:00:13
estimated time: 00:00:35
-----
progress: 30%
elapsed time: 00:00:16
estimated time: 00:00:37
-----
progress: 30%
elapsed time: 00:00:16
estimated time: 00:00:37
-----
progress: 32%
elapsed time: 00:00:17
estimated time: 00:00:36
-----
progress: 32%
elapsed time: 00:00:22
estimated time: 00:00:46
-----

```

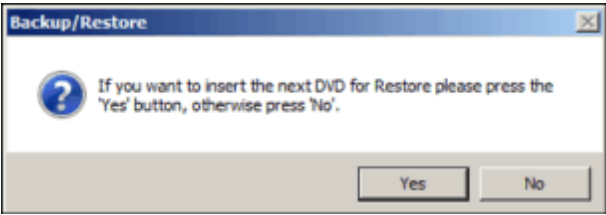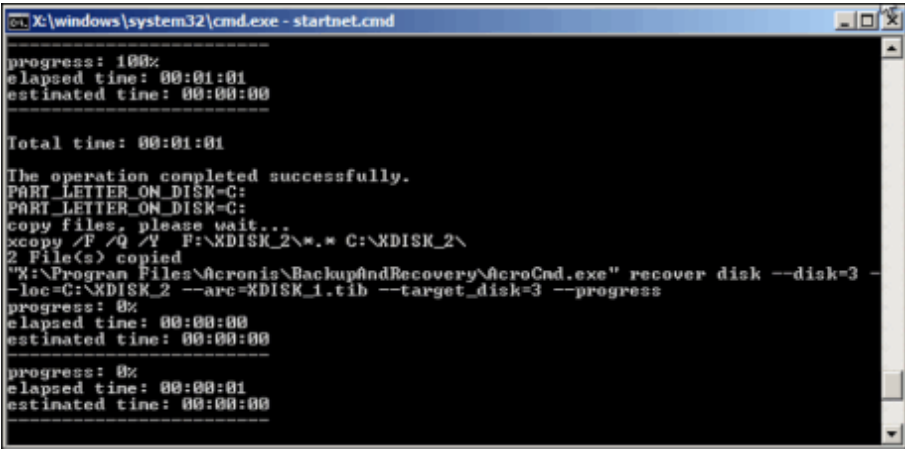
ACTION	MRAWP
Wait until the following popup appears. Fig. 82: Next DVD request popup 	☐
Confirm the above shown "Next DVD request popup" with "No".	☐
Now the last disk is restored from the image. Fig. 83: Last disk restore 	☐
When all data have been successfully restored, the following popup is displayed. Fig. 84: Backup/Restore successfully finished popup 	☐
Remove the ejected medium, close the drive, and confirm this popup with OK.	☐
The system is rebooted.	☐
After the reboot, an autologin as administrator is executed.	☐
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 10 minutes. Note: Don't accept any displayed popup or press a button on any displayed popup.	☐

ACTION

MRAWP


During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:

Fig. 85: Displayed popup

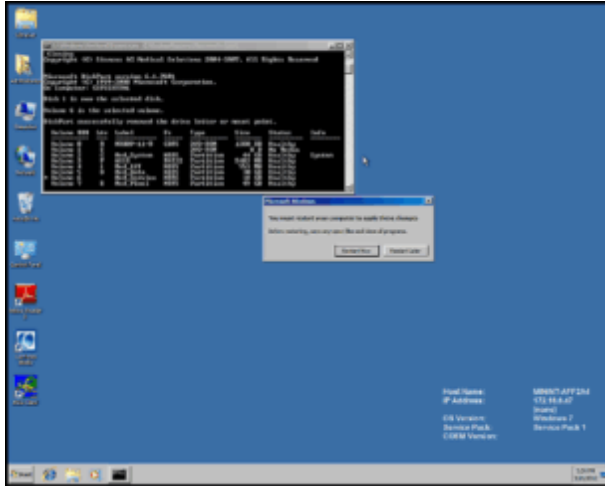


Note: Don't accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don't touch or accept it. Keep it displayed on the screen. It will be closed automatically.


Fig. 86: Windows popup

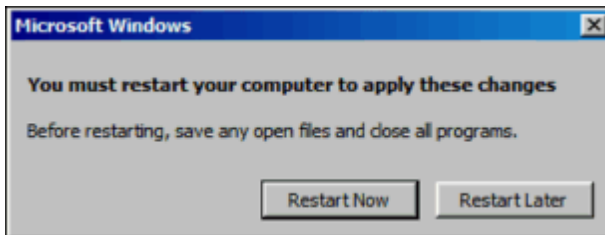


Note: The following popup can be ignored too, if it is displayed.

Fig. 87: Network location popup

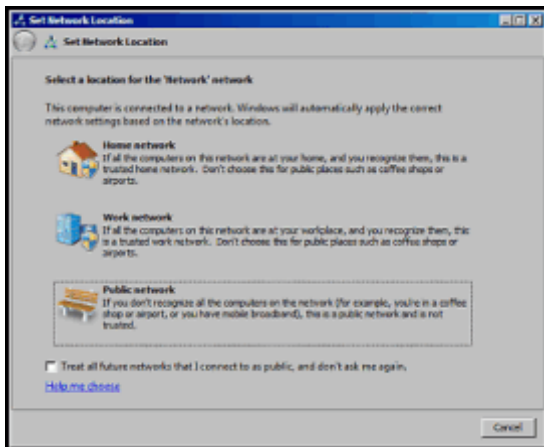

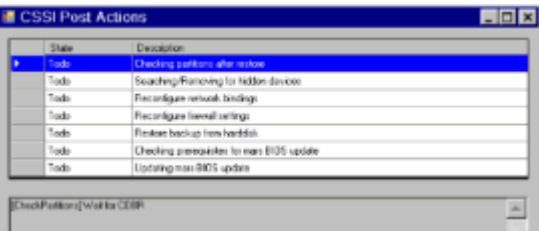ACTION		MRAWP
The CSSI Post popup is displayed (see (→ Fig. 88 Page 76)) Fig. 88: CSSI Post Action popup 		
The following CSSI Post Action steps are performed: ■ Checking partitions after restore (duration about 2 min.) ■ Searching/Removing for hidden devices (duration about 1 min.) ■ Reconfigure network bindings (duration about 3 min.)		
■ The system is rebooted automatically. No user interaction is necessary.		☐
■ A autologin as user "administrator" is performed and the CSSI Post Actions popup (→ Fig. 88 Page 76) is displayed. The following CSSI Post Action steps are performed: - Reconfigure firewall settings (duration about 1 min.) - Restore backup from hardisk (only IP settings, duration about 1 min.)		☐
■ The system is rebooted automatically. No user interaction is necessary.		☐

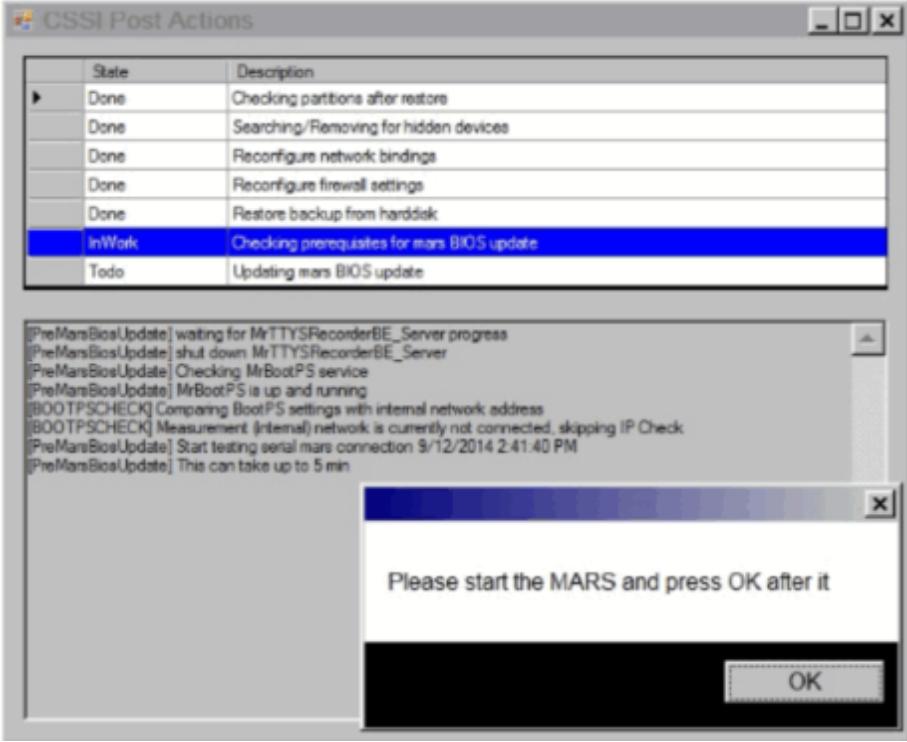ACTION	MRAWP
■ A autologin as user "administrator" is performed and the CSSI Post Actions popup (→ Fig. 88 Page 76) is displayed. The following CSSI Post Action steps are performed: - Checking prerequisites for mars BIOS update Wait until you get the Please start the MARS... popup (duration about 5 min.). <i>Fig. 89: MARS reboot pop up</i> 	☐
■ Connect the power line connector from the MARS in the Electronic Power Cabinet.	☐
■ Confirm the Please start the MARS... popup with OK.	☐
■ The MARS is booted up and the following CSSI post action step continues: - Checking prerequisites for mars BIOS update (duration about 15 min.) Automatically a reboot is performed after this step is done.	☐

ACTION	MRAWP
Note: Is the reboot not performed automatically, please perform the following workaround. Run the CSSI Post Actions a second time. - Close the CSSI Post Actions by pressing the X. - Start a cmd tool window (Start ⇒ All Programs ⇒ Command Prompt). - In the cmd tool window change the directory to <code>C:\MedCom\utils</code> - Execute again the CSSI Post actions by typing in the command: <code>CSSIPostRestoreProcess.exe ↵</code> Now the CSSI Post actions starts again. This takes about 20 minutes. After this time a reboot is done automatically. If the CSSI Post actions tool fails again, please contact the HSC MR in Erlangen.	
The system performs a autologin as meduser..	☐



In order to start the CSSI Post actions manual (in case of problems) do the following:
 Start a cmd tool window (Start ⇒ All Programs ⇒ Command Prompt). In the cmd tool window change the directory to `C:\MedCom\utils`
 and execute again the CSSI Post actions by typing in the command: `CSSIPostRestoreProcess.exe ↵`

3.1.1.8 MRAWP (MRC) System Start-up



If the customer has TIM docking table (mobile), the table must be connected during the whole firmware/loadware update process.

Is the TIM docking table removed during the firmware/loadware update process the PSY (Power Stage Y) is destroyed.

ACTION	MRAWP
Wait until the main UI is displayed.	☐


ACTION

MRAWP

- The firmware update starts, and a the **MR Scanner & Acquisition** message popup is automatically displayed.




Fig. 90: Example for MR Scanner & Acquisition popup

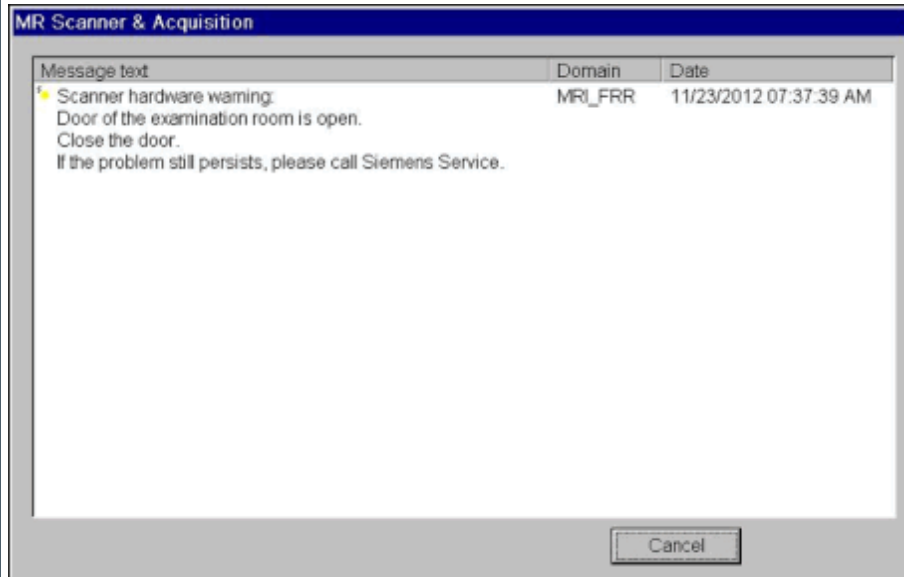

Note: A message is displayed for each component when the download is started. Also the duration of the download is given in this message. Below is an example of it:

Fig. 91: Message popup

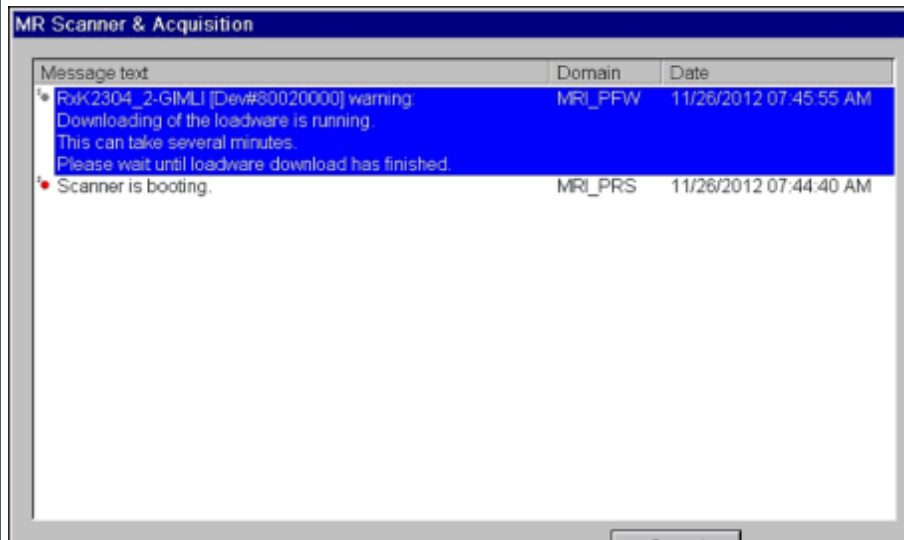

This message will disappear automatically when the component download is finished. -> You can check if recently an update is running.

Note: For the PTAB an error is displayed until the firmware was loaded successful.

If no message is displayed the wrong language setting is used. Please correct it according (→ Correct syngoUI Language Settings / Page 180) after the firmware update. In this case please observe the firmware update in the windows eventlog.

ACTION	MRAWP												
The following popup be displayed as an information, that the MARS will be rebooted automatically, Press OK to close the popup. Fig. 92: Reboot message popup MRScanner Auto rebooting scanner after loadware download. This is auto reboot try number 1. Please wait approx. 10 min until boot has finished. OK													
Wait until the component firmware update has finished. The firmware update has finished when the following message is displayed in the MR Scanner & Acquisition message popup: - Feedback Loop Error: Feedback Loop Parameters do not match the current TX-Box hardware configuration. Please restart Meas & Recon or switch off/on the Scanner. If the problem persists, please call Siemens Service. The reason for this error message is the new Tx-Box firmware which resets the parameters. This can take up to 60 minutes max. until you get this message! Fig. 93: Message popup "Feedback Loop Error" MR Scanner & Acquisition <table><tr><th>Message text</th><th>Domain</th><th>Date</th></tr><tr><td>Scanner is not ready for measurement. Reboot Scanner, switch off/on Scanner. Reboot syngo MR. Please save System Log-Files immediately by selecting "System/Control/Tools" and "Save System Log-Files". Please call Siemens Service.</td><td>MRI_PRS</td><td>11/26/2012 08:15:26 AM</td></tr><tr><td>Booting the Scanner has failed. See previous messages for details.</td><td>MRI_PRS</td><td>11/26/2012 08:15:26 AM</td></tr><tr><td>Feedback Loop Error: Feedback Loop Parameters do not match current Tx Box hardware configuration. Please restart Meas & Recon or switch off/on Scanner. If the problem still persists, please call Siemens Service.</td><td>MRI_FRR</td><td>11/26/2012 08:15:26 AM</td></tr></table> Cancel	Message text	Domain	Date	Scanner is not ready for measurement. Reboot Scanner, switch off/on Scanner. Reboot syngo MR. Please save System Log-Files immediately by selecting "System/Control/Tools" and "Save System Log-Files". Please call Siemens Service.	MRI_PRS	11/26/2012 08:15:26 AM	Booting the Scanner has failed. See previous messages for details.	MRI_PRS	11/26/2012 08:15:26 AM	Feedback Loop Error: Feedback Loop Parameters do not match current Tx Box hardware configuration. Please restart Meas & Recon or switch off/on Scanner. If the problem still persists, please call Siemens Service.	MRI_FRR	11/26/2012 08:15:26 AM	
Message text	Domain	Date											
Scanner is not ready for measurement. Reboot Scanner, switch off/on Scanner. Reboot syngo MR. Please save System Log-Files immediately by selecting "System/Control/Tools" and "Save System Log-Files". Please call Siemens Service.	MRI_PRS	11/26/2012 08:15:26 AM											
Booting the Scanner has failed. See previous messages for details.	MRI_PRS	11/26/2012 08:15:26 AM											
Feedback Loop Error: Feedback Loop Parameters do not match current Tx Box hardware configuration. Please restart Meas & Recon or switch off/on Scanner. If the problem still persists, please call Siemens Service.	MRI_FRR	11/26/2012 08:15:26 AM											

3.1.1.9 MRAWP (MRC) Tune up and QA Steps

ACTION	MRAWP
Restore the Tuneup status: Go in local service UI: - Select Utilites ⇒ Escape to OS and type in the command <code>MrSesoRestoreTuneUpStatus.bat</code> and press <ENTER> Fig. 94: SESO output  This step is restoring the last TuneUp&QA steps. The QA steps for the surface coils are restored too.	☐
Click Home and select Tune Up.	☐
Perform the following tune-up step first (keep the default selection of the tune up setps): - Feedback Loop Calibration During this step the MARS is rebooted.	☐
Perform all other tune-up steps that are set to the "ToDo" status	☐
Note: If you have a satellite console MRWP (MRSC), you can start with the step (→ MRWP (MRSC) Preliminary Steps before the Software Installation / Page 93) now in parallel to the tune up and QA.	
Perform all QA steps that are set to the "ToDo" status	☐

3.1.1.10 MRAWP (MRC) Surface Coils QA Steps



For all coils used before the software update, it is not necessary to perform a Coil QA. The QA status is restored.

3.1.1.11 MRAWP (MRC) Network Printer Installation

ACTION	MRAWP
Install the network paper printer according to the software installation manual (this step is not part of the dynamic data) if previously installed on the MRAWP or MRWP.	☐
If applicable, set the installed printer as the default printer. Select Start -> Settings -> Printers and Faxes Select the printer and click the right mouse button. In the menu displayed, select the "Default printer" entry.	☐

3.1.1.12 MRAWP (MRC) Function Test of the DICOM Nodes

ACTION	MRAWP
Send an image from the Patient Browser to a all connected DICOM nodes. The images have to be sent successfully	☐
Send a film sheet to all connected DICOM cameras. The film sheet has to be sent successfully to the DICOM camera.	☐
Retrieve a test patient from the RIS/HIS system, if available. The patient should be available on the system.	☐



The new DICOM Node Configuration Option "Supports Presentation State" must not be enabled for any DICOM storage node. At the moment the NUMARIS/4 VE11 does not support the "Presentation State" for created MR images.

3.1.1.13 MRAWP (MRC) System Management



Due to the CSSI process, system management cannot be activated at the factory. System management has to be activated at the customer site.

Installation and Activation Steps

ACTION	MRAWP
Select Options --> Service -> Local Service.	☐

ACTION	MRAWP
Enter the valid service key.	☐
In the local Service UI, select Configuration	☐
Check that the settings for service region (Europe, USA, or Asia) under Local Host Site Info are correct. Correct the settings, if necessary, and save the changes. This entry assigns the IP addresses for the RDIAG server and the System Management server.	☐
Under the "Local Host" Routing, type in the IP address of the router or gateway used for SRS. In some cases, this entry is the default gateway that is used for the TCP/IP LAN configuration. See logfile C:\MedCom\MriSiteData\ExtBaRe\routing.txt for the routing entry. For the installation of System Management, the correct entry under routing is essential. If this routing entry is not set (empty), you will get an error popup when entering the System Management configuration.	☐
Select System Options. Fig. 95: Configuration 	☐


ACTION	MRAWP
■ Under "List of System Options", click the System Management checkbox Fig. 96: System options The screenshot shows a web browser window titled "Siemens Med Service Software - Windows Internet Explorer". The address bar shows "Configuration". The page has a navigation bar with "Service Key Expiry Date: 14-May-2012", "Reports", "Home", and "Help". Below the navigation bar, there is a checkbox labeled "enable help pages". The main content area is titled "System type" with a dropdown menu showing "Skysa". Below this is a section titled "List of system options" with a list of checkboxes: "Attached to Network" (checked), "Internet Based Connectivity" (unchecked), "DICOM Print Devices" (checked), "DICOM Offline Devices" (checked), "DICOM Networking" (checked), "DICOM HIS/RIS" (checked), "Image Import/Export" (checked), "Paper Printer" (unchecked), "EPRI" (unchecked), "System Management" (checked), "MR MBOX" (unchecked), and "CAN24" (unchecked). At the bottom of the window, there is a "Next" button.	☐
■ Select Next.	☐

ACTION

MRAWP

- On the System Management Master configuration page


Fig. 97: System Management master page

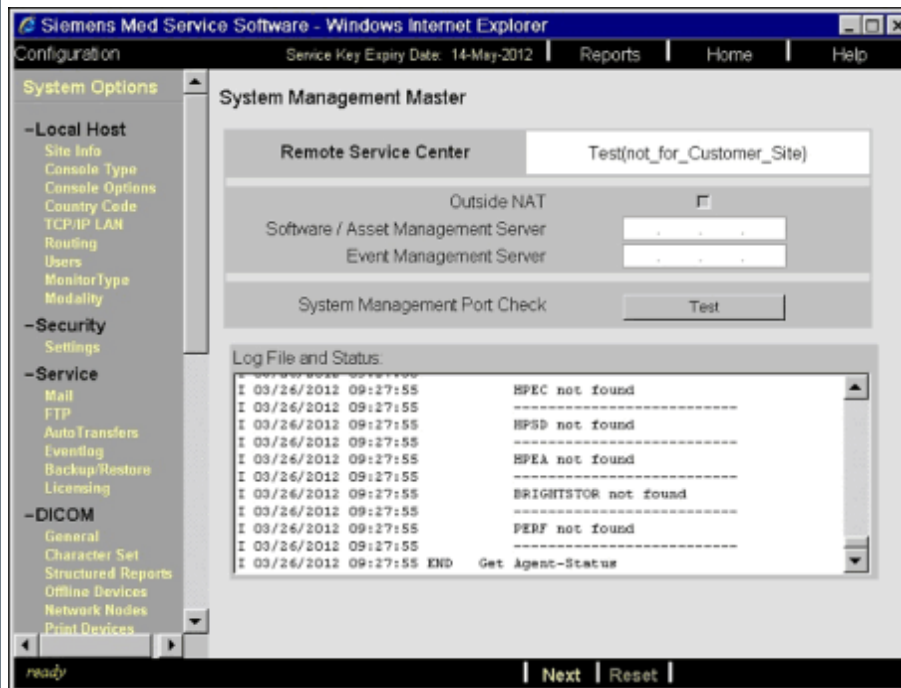

See the report file

C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html
for the Outside NAT configuration.

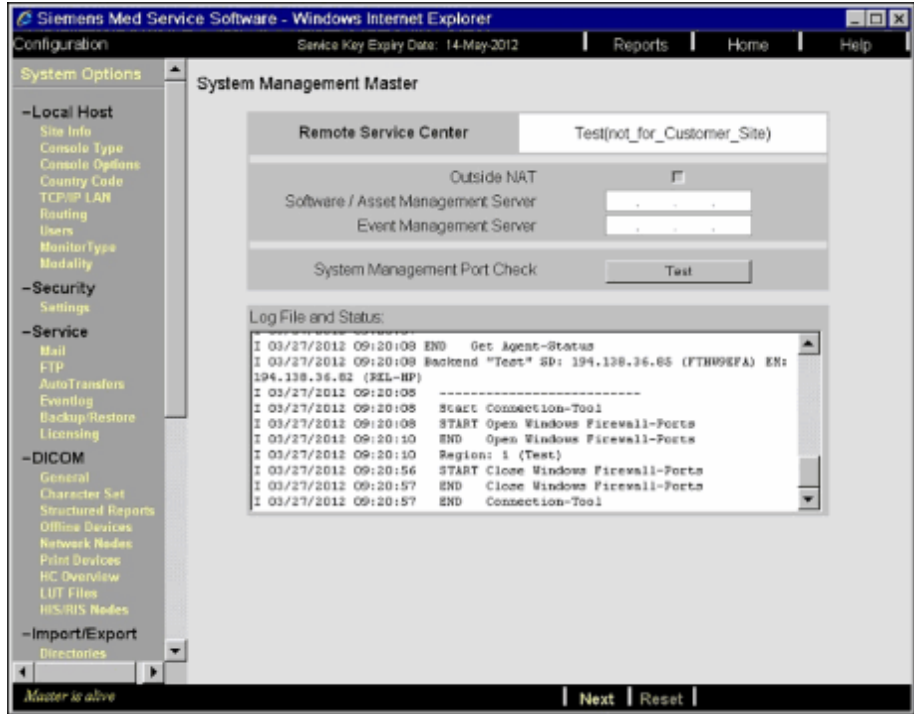

ACTION	MRAWP
■ Use the Test button to check whether the System Management Server can be accessed. Fig. 98: System Management master can be reached  The log file for the test can be found under the directory: C:\MedCom\service\mwtools\CCT\Logs	☐
In the event of problems with the connection test, contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	
■ Select Next.	☐


ACTION

MRAWP

- Under **System Mgmt. Agent Controls**, set the following configuration:
Mark the checkboxes "Go Active" and "Go Online".
(The other checkboxes have already been marked).

☐

Fig. 99: System Management activation

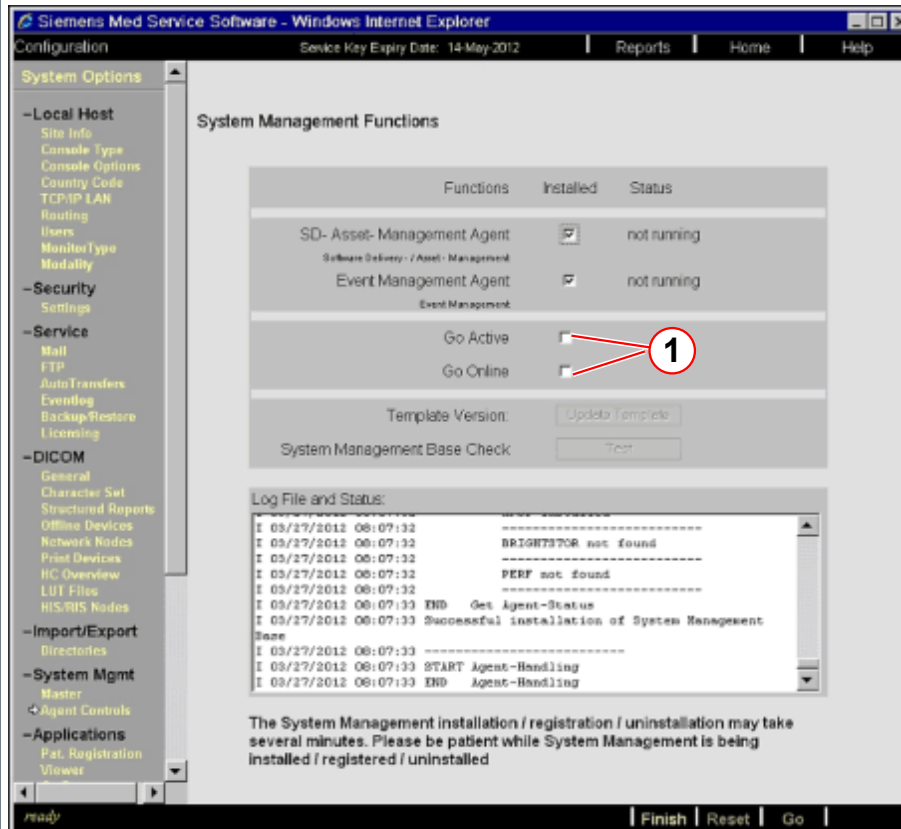

(1) Mark both checkboxes

In the Log File and Status area, the current status of System Management activation/installation is displayed.

- Click **Go**.

☐

Note: The installation/registration will take approximately 10 to 15 minutes.

ACTION

MRAWP

- If activation is successful, the following is displayed:

Fig. 100: System Management function activation OK


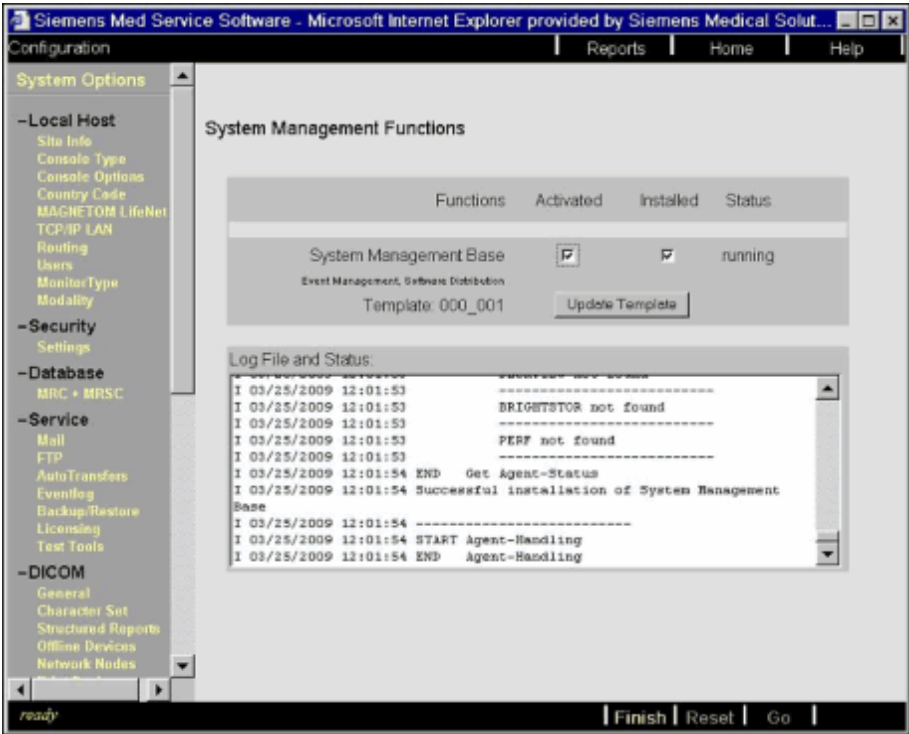


Fig. 101: System Management activation OK

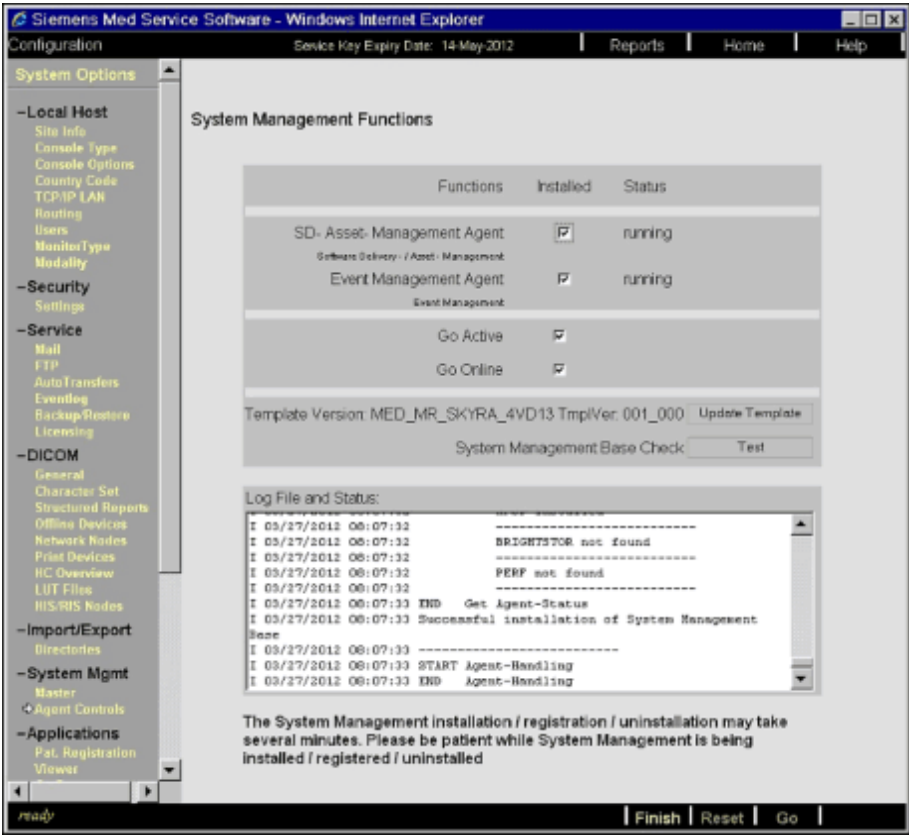

Select **Finish** to continue

ACTION	MRAWP
Click Finish and close the local service UI.	☐
In the event of problems (installation/activation fails a second time), contact the SRS help desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	


3.1.1.14 MRAWP (MRC) Restoring changed Customer Measurement Programs from VD13A/VD13B/VD13C/VD13D/VD13E



The following step restores all changed customer protocols. Only these protocols will be converted after the reboot.

ACTION	MRAWP
Insert the MasterBackup CD labeled Master Backup VD13A or VD13B or VD13C or VD13D or VD13E MRAWP (actual date) created see (→ Performing a Master Backup on the MRC, for changed Customer Protocols / Page 60) into the CD/DVD-ROM drive R:	☐
Select Options --> Service --> Local Service.	☐
Enter the valid service key and click OK.	☐
Select Backup&Restore in the service platform.	☐
Select under Command Restore	☐
Select under Drives [-R-] Is the drive letter not correct please perform (→ Drive Letter Assignment for the CD/DVD ROM Drive and/or for the CD/DVD-R Drive / Page 178).	☐
Select under Archive the MasterBackup and under Groups the entry N4Custom-er.	☐
Press Go to start the restore. Wait until the restore has finished.	☐
Remove the MasterBackup CD labeled Master Backup VD13A or VD13B or VD13C or VD13D or VD13E MRAWP (actual date) from the CD/DVD-ROM drive R:	☐
Click Home and confirm the requested reboot message popup with OK.	☐

3.1.1.15 MRAWP (MRC) Performing a Win 7 Image Backup on the MRAWP (MRC)

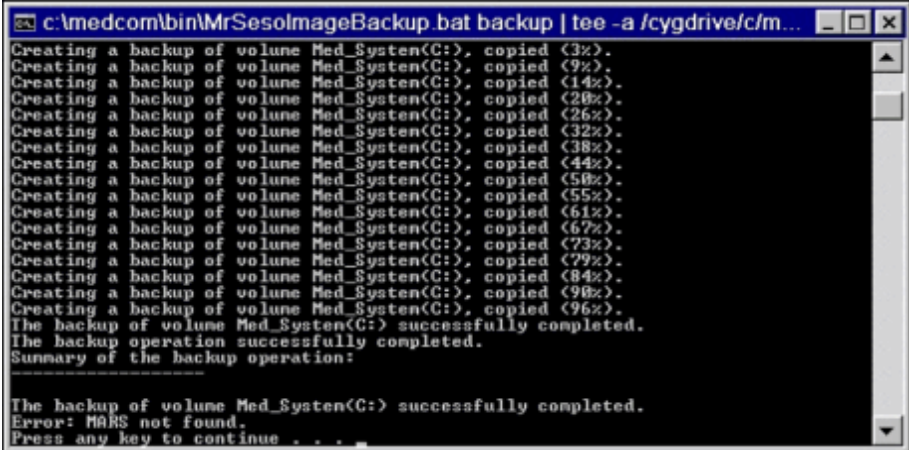ACTION	MRAWP
■ Create the Win 7 image backup of the system disk - Enable advanced user. - Start -> All Programs - Select Image Backup The administrator password is requested first: <i>Fig. 104: Password request popup</i>  Type in the administrator password (the input is not displayed) and press ↵ Now the Windows 7 image is automatically created from the system disk, and a copy is transferred to MARS. Wait until the procedure is finished. This takes about 10 minutes. <i>Fig. 105: Win7 image Backup</i>  Close the Image Backup window by pressing <SPACE>	☐

3.1.1.16 MRAWP (MRC) Enable the Auto Login Function



The following steps are required only if auto login is not enabled.

ACTION	MRAWP
■ Select Options --> Service --> Local Service .	☐
■ Enter the valid service key and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select local host Users	☐
■ Mark the checkbox " Enable Auto login " to enable the auto login function.	☐
■ Click Save .	☐
■ Click Home to leave the configuration.	☐
■ Close the local service UI.	☐

3.1.1.17 MRAWP (MRC) Additional steps for TIM docking table (mobile)



The following steps are necessary, if the customer has more than one TIM docking table (mobile). The new firmware/loadware has to be loaded before the table can be used in the new Numaris version VE11A or VE11B.

With the new firmware/loadware the docking table can not be used on MR Systems with former software version (VD13A/VD13B/VD13C/VD13D/VD13E).

ACTION	MRAWP
■ Connect the other docking table to the MR system	☐
■ A reboot of the MARS is required to start the update of the docking table. Perform the reboot as follows: - In the syngo main UI, select System - Control. - In the System Manager popup, select the Meas & Recon tab. - Press the Reboot button under Meas & Recon Control. - Confirm the displayed popup with Yes. MARS is now booting.	☐
■ Wait until the docking table firmware update has finished. Observe the firmware/loadware update in the MR Scanner & Acquisition message popup. The update is finished when the Measurement system is ready for measurement.	☐




Is no MRWP (MRSC) available the following steps can be skipped. Go ahead and continue with (→ Final Check - to be performed for MRC and MRSC (if applicable) / Page 110)

3.1.1.18 MRWP (MRSC) Preliminary Steps before the Software Installation

ACTION	MRWP
- Be sure that the network cables at the back of the console are connected. MRSC: One network cable at the back of the computer.	☐
Be sure that any device (e.g., monitor, in-room MRC, etc.) connected to the second video output of the host is disconnected. Disconnect any additionally connected USB device (keep only the USB mouse and the USB keyboard connected).	☐
If applicable, check whether the local connected printer (connected to the parallel port) is switched ON. During the installation, the printer is automatically detected and the correct printer driver is installed.	☐
Switch on the MRSC.	☐
Log in as user "administrator".	☐



3.1.1.19 MRWP (MRSC) BIOS Update on MRWP (MRSC)




ACTION	MRWP
Insert the first DVD of the delivered CSS image set labeled COBR Recovery Image MR system software VE11A or VE11B MR system serial no: <x> MR WP into the CD/DVD burner drive (S:\).	☐
Insert the delivered CD labeled with "MR-VD BIOS CD VE11.. " into the host CD/DVD-ROM Drive" (R:\)	☐
Invoke a command shell from the start menu (Start ⇒ All Programs⇒ Accessories ⇒ Command Prompt)	☐
In the command shell, type in: S:\Start.bat	☐



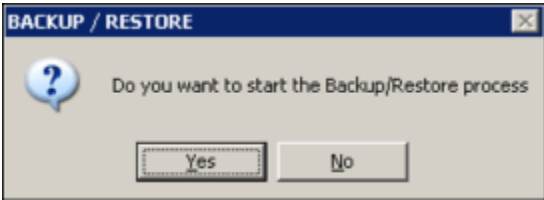
ACTION	MRWP
Confirm the Installation popup with "Yes" Fig. 106: Backup/Restore popup 	☐
The running script performs a reboot at the end. This may take some time. Please wait!	☐
The BIOS CD will be booted automatically.	☐
The BIOS update runs automatically, no user action is necessary, and ends when the following message is displayed: <pre>BIOS Installation Done - System Ready for Software Installation Please remove BIOS CD ! Please Turn Off Power or Reset System ! C:\></pre>	☐
Remove the delivered BIOS CD from the CD/DVD drive (R:\)	☐
Remove the first DVD of the delivered CSS image set labeled COBR Recovery Image MR system software VE11A or VE11B MR system serial no: <x> MR WP from the CD/DVD burner drive (S:\).	☐
Insert the first DVD of the delivered CSS image set labeled COBR Recovery Image MR system software VE11A or VE11N MR system serial no: <x> MR WP into the CD/DVD-ROM drive (R:\).	☐
Press key combination <CTRL> + <ALT> + (on a German keyboard use <Strg> + <Alt> + <Entf>) to reboot and continue with "CSSI installation MRC".	☐




3.1.1.20 MRWP (MRSC) CSSI Installation

ACTION	MRSC
A reboot is performed and the COBR tool boots.	☐

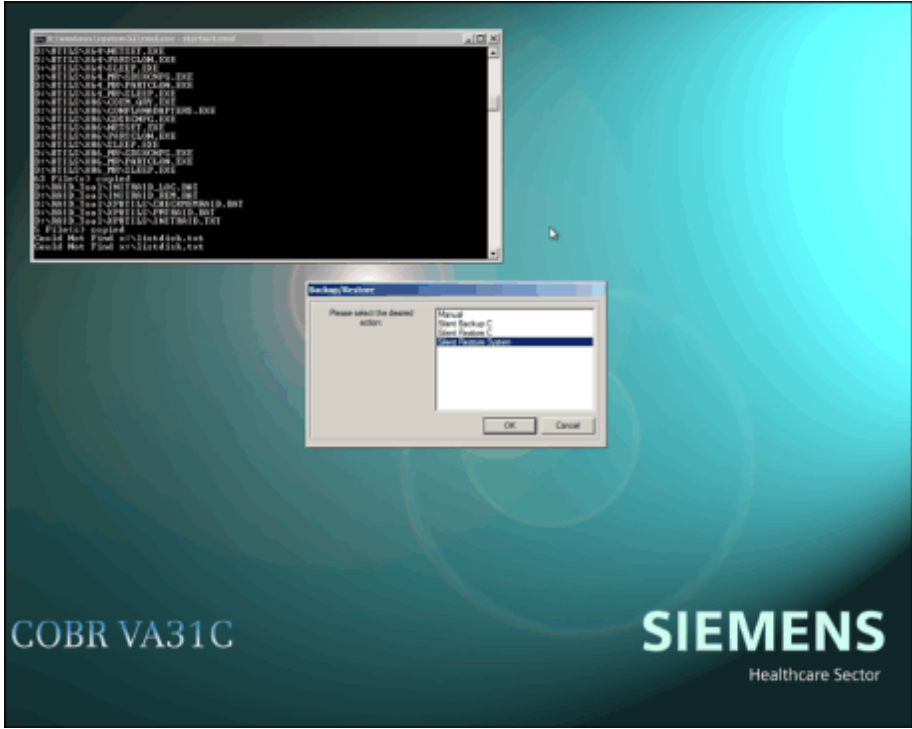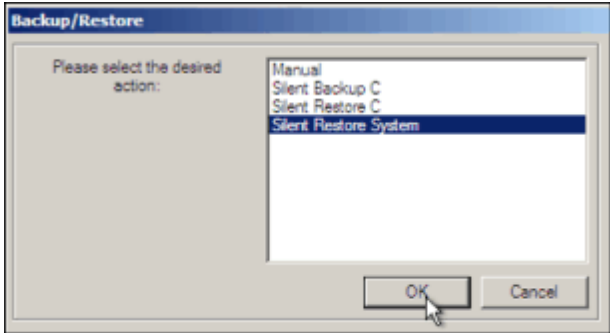ACTION	MRSC
Note: If you get the following error message "Operating system not found", please remove the first DVD of the delivered CSS image set from the CD/DVD drive (R:) and insert it into the CD/DVD burner drive (S:) and perform a reboot.	
Wait until the COBR tool window is displayed. Fig. 107: COBR tool window 	☐
In the Backup/Restore popup, select <code>Silent Restore System</code> and press OK Fig. 108: Backup/Restore popup Silent Restore System 	☐

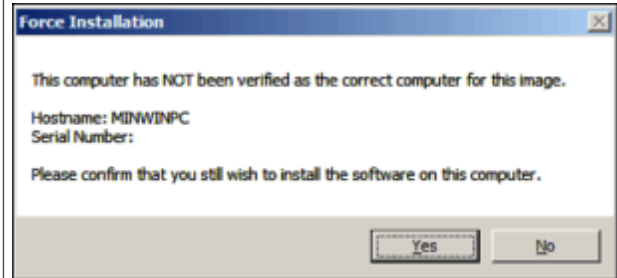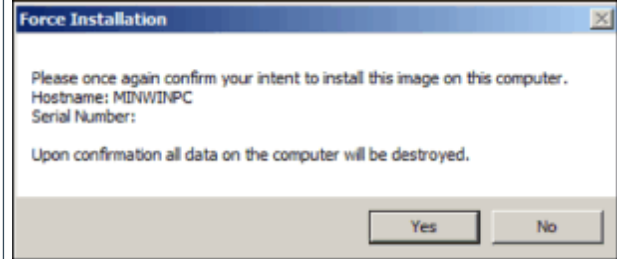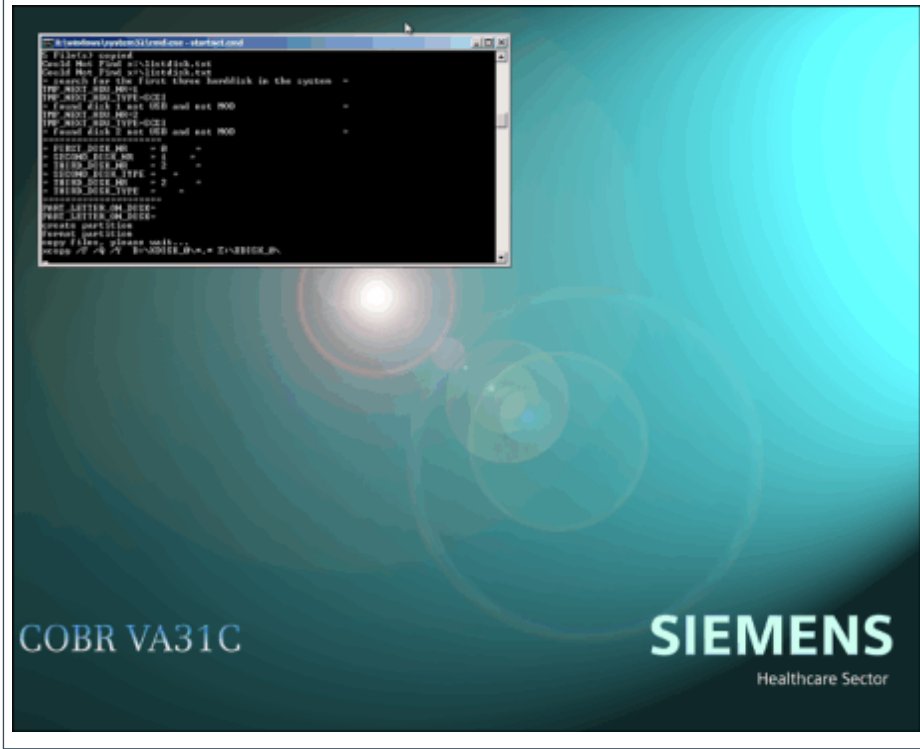ACTION	MRSC
■ The following popup is displayed (example is shown below) <i>Fig. 109: Force Installation popup</i> A Windows-style dialog box titled "Force Installation". The text inside reads: "This computer has NOT been verified as the correct computer for this image." followed by "Hostname: MINWINPC" and "Serial Number:". Below this, it says "Please confirm that you still wish to install the software on this computer." At the bottom are "Yes" and "No" buttons. Confirm this popup with Yes.	☐
■ Confirm the Force Installation popup (example is shown below) with Yes. <i>Fig. 110: Force Installation popup</i> A Windows-style dialog box titled "Force Installation". The text inside reads: "Please once again confirm your intent to install this image on this computer." followed by "Hostname: MINWINPC" and "Serial Number:". Below this, it says "Upon confirmation all data on the computer will be destroyed." At the bottom are "Yes" and "No" buttons. This Force Installation popup is displayed twice. Confirm it each time with Yes.	☐
■ The system restore process starts now. The following shows an example of the screen display. <i>Fig. 111: COBR restore process start</i> The screen shows a Siemens COBR VA31C restore process. The background is a blue gradient with a Siemens logo and "Healthcare Sector" text. In the top left, a command prompt window displays the output of the "start.bat" script, showing file copying and disk formatting steps. The text "COBR VA31C" is visible in the bottom left corner.	☐

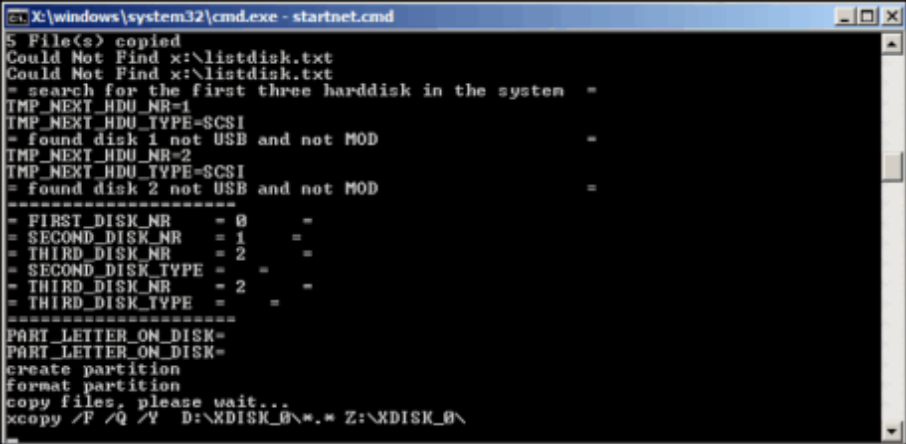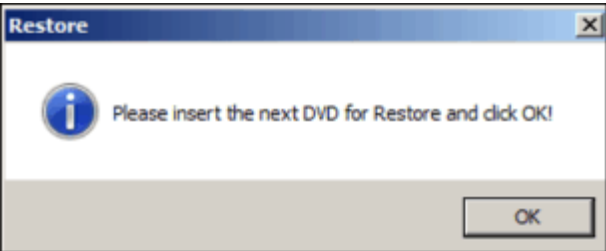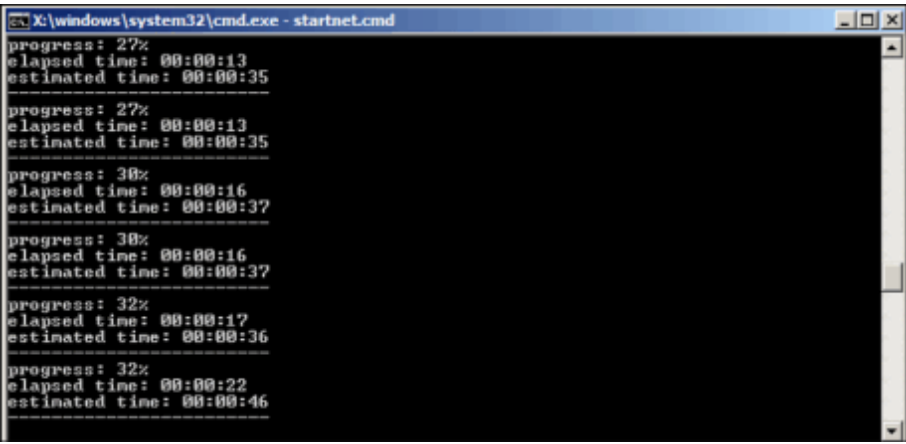ACTION	MRSC
First all image files from the DVD for the system disk (disk1) are copied from the DVD to the hard disk. Fig. 112: Copy image files  After the files are copied, the DVD is automatically ejected.	☐
Remove the ejected medium from the CD/DVD-ROM drive, insert the next DVD from the system image DVD set, and close the drive.	☐
Confirm the request popup with OK. Fig. 113: Next DVD request popup 	☐
The copy procedure continues. The next DVD is requested automatically until all required data are copied. The copy and the confirmation of the popup (→ Fig. 113 Page 97) have to be repeated for all available CSSI DVDs.	
Now the image restore starts. Fig. 114: Image restore progress output 	☐

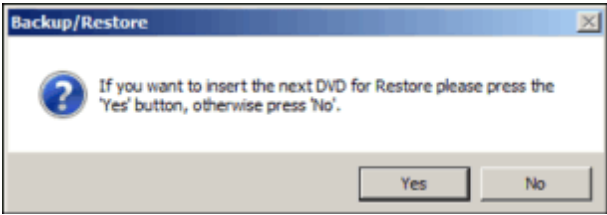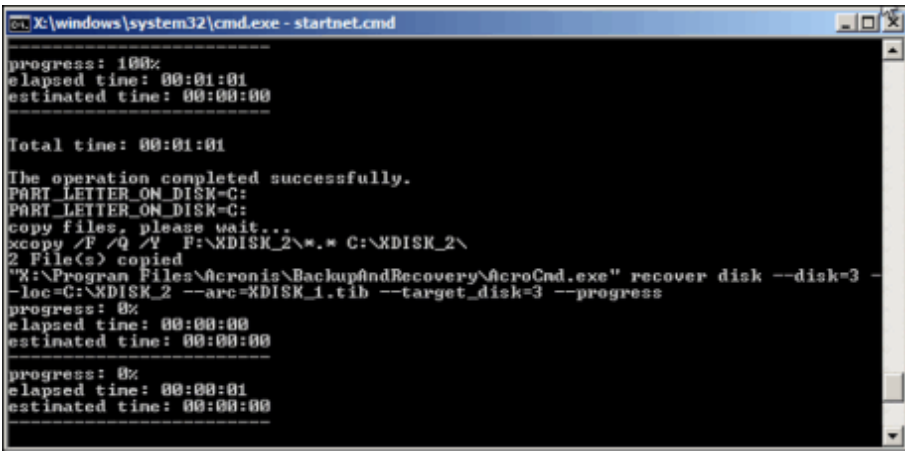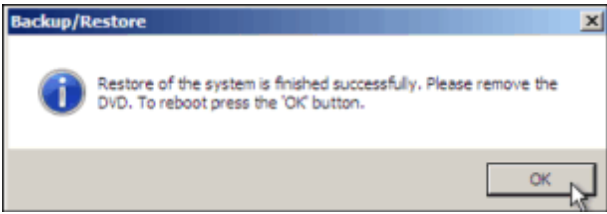
ACTION	MRSC
Wait until the following popup appears. Fig. 115: Next DVD request popup 	☐
Confirm the above shown "Next DVD request popup" with "No".	☐
Now the last disk is restored from the image. Fig. 116: Last disk restore 	☐
When all data have been successfully restored, the following popup is displayed. Fig. 117: Backup/Restore successfully finished popup 	☐
Remove the ejected medium, close the drive, and confirm this popup with OK.	☐
The system is rebooted.	☐
After the reboot, an autologin as administrator is executed.	☐
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 5 minutes. Note: Don't close or confirm any displayed popups.	☐

ACTION

MRSC

During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:


Fig. 118: Displayed popup

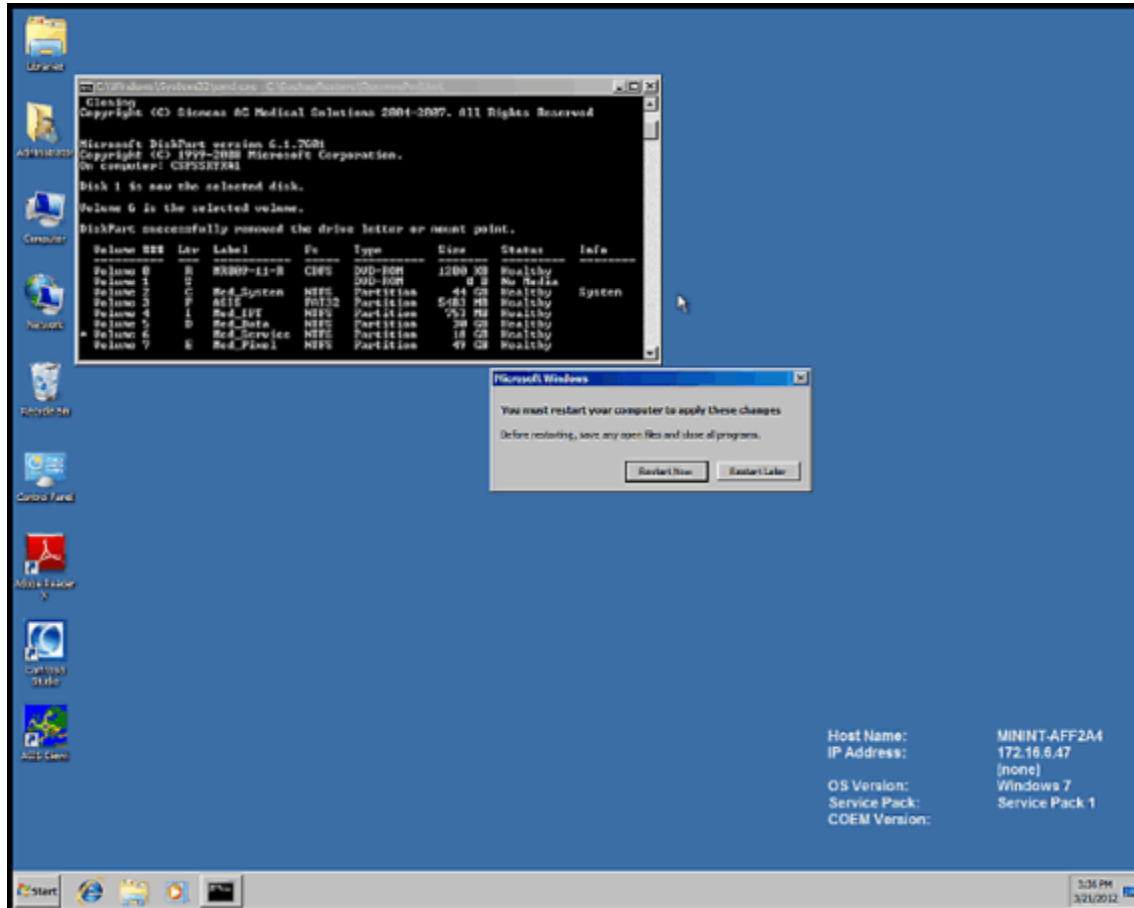



The CSSI Post popup is also displayed (see Fig. 60 “CSSI Post Actions popup”)

Note: Don’t accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don’t touch or accept it. Keep it displayed on the screen. It will be closed automatically.

Fig. 119: Windows popup

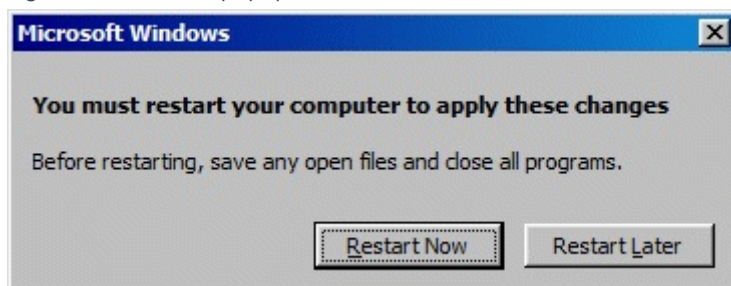

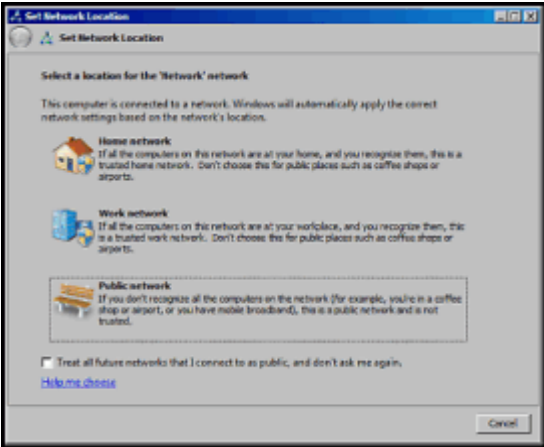ACTION		MRSC
Note: The following popup can be ignored too, if it is displayed. Fig. 120: Network location popup 		
■ The system is rebooted.		☐
■ A autologin as user "administrator" is performed and the following CSSI post popup is displayed.		☐

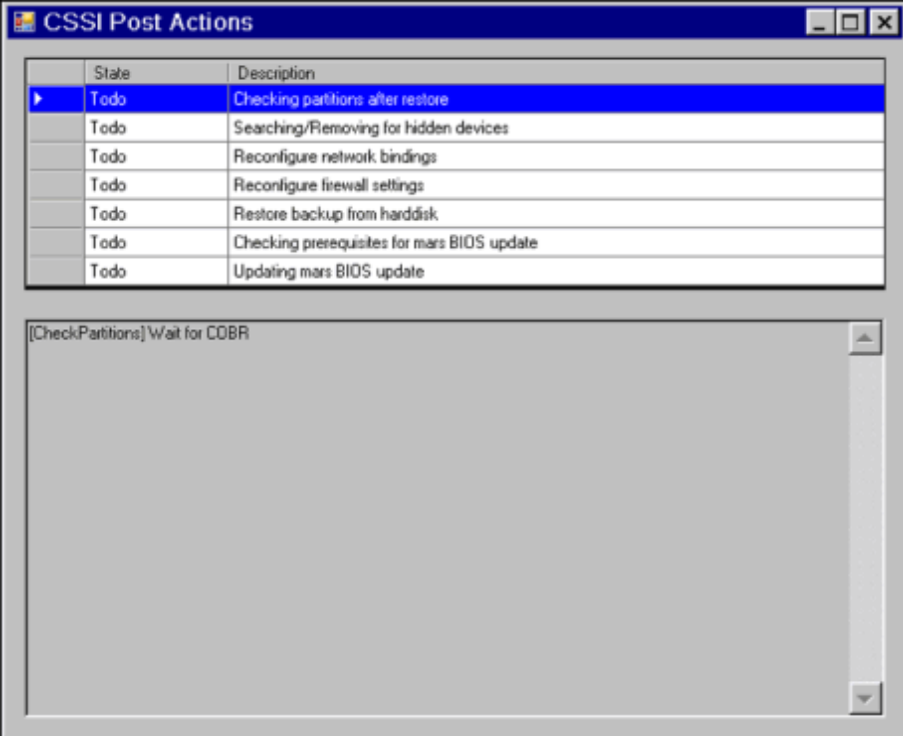ACTION	MRSC
The CSSI post popup shows the status and progress of the different post CSSI update steps. <i>Fig. 121: CSSI Post Actions popup</i>  Please wait until the steps are finished. During these post CSSI update steps the host will be rebooted automatically. Note: No user interaction is necessary.	☐
The system performs a autologin as meduser after all post CSSI update steps are done.	☐

3.1.1.21 MRWP (MRSC) Checking the Patient Browser Connection with Shared Database

ACTION	MRWP
Use the Patient Browser to check whether the MRC database can be reached.	☐

3.1.1.22 MRWP (MRSC) Network Printer Installation on the MRSC

ACTION	MRWP
Install the network paper printer according to the software installation manual (this step is not part of the dynamic data) if previously installed on the MRAWP or MRWP.	☐
If applicable, set the installed printer as the default printer. Select Start -> Settings -> Printers and Faxes Select the printer and click the right mouse button. In the menu displayed, select the "Default printer" entry.	☐

3.1.1.23 MRWP (MRSC) Function Test of the DICOM Nodes

ACTION	MRWP
Send an image from the Patient Browser to all connected DICOM nodes. The images have to be sent successfully	☐
Send a film sheet to all connected DICOM cameras. The film sheet has to be sent successfully to the DICOM camera.	☐



The new DICOM Node Configuration Option "Supports Presentation State" must not be enabled for any DICOM storage node. At the moment the NUMARIS/4 VE11 does not support the "Presentation State" for created MR images.


3.1.1.24 MRWP (MRSC) System Management



Due to the CSSI process, system management cannot be activated at the factory. System management has to be activated at the customer site.

Installation and Activation Steps

ACTION	MRWP
Select Options --> Service -> Local Service.	☐
Enter the valid service key.	☐
In the local Service UI, select Configuration	☐
Check that the settings for service region (Europe, USA, or Asia) under Local Host Site Info are correct. Correct the settings, if necessary, and save the changes. This entry assigns the IP addresses for the RDIAG server and the System Management server.	☐

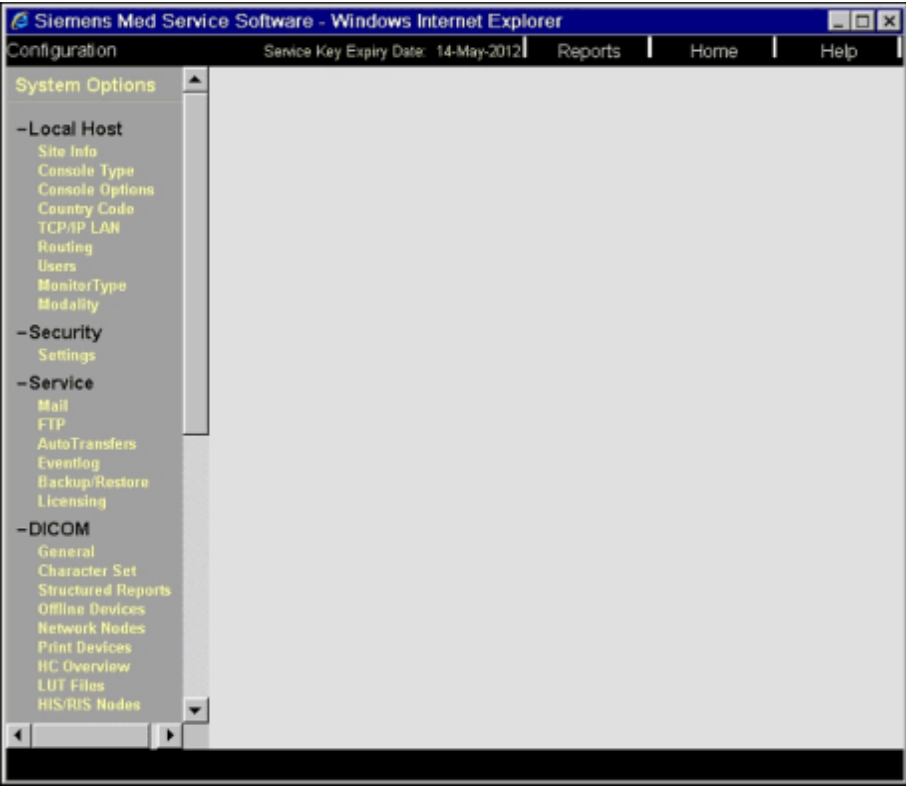ACTION	MRWP
Under the "Local Host" Routing, type in the IP address of the router or gateway used for SRS. In some cases, this entry is the default gateway that is used for the TCP/IP LAN configuration. See logfile C:\MedCom\MriSiteData\ExtBaRe\routing.txt for the routing entry. For the installation of System Management, the correct entry under routing is essential. If this routing entry is not set (empty), you will get an error popup when entering the System Management configuration.	☐
Select System Options. Fig. 122: Configuration 	☐

ACTION	MRWP
■ Under "List of System Options", click the System Management checkbox Fig. 123: System options The screenshot shows a web browser window titled "Siemens Med Service Software - Windows Internet Explorer". The address bar shows "Configuration". The page has a navigation bar with "Service Key Expiry Date: 14-May-2012", "Reports", "Home", and "Help". Below the navigation bar, there is a checkbox labeled "enable help pages". The main content area is titled "System type" with a dropdown menu showing "Skyra". Below this is a section titled "List of system options" containing a list of options with checkboxes: "Attached to Network" (checked), "Internet Based Connectivity" (unchecked), "DICOM Print Devices" (checked), "DICOM Offline Devices" (checked), "DICOM Networking" (checked), "DICOM HIS/RIS" (checked), "Image Import/Export" (checked), "Paper Printer" (unchecked), "EPRI" (unchecked), "System Management" (checked), "MR MBOX" (unchecked), and "CAN24" (unchecked). At the bottom of the window, there is a "Next" button.	☐
■ Select Next.	☐


ACTION

MRWP

- On the System Management Master configuration page




Fig. 124: System Management master page

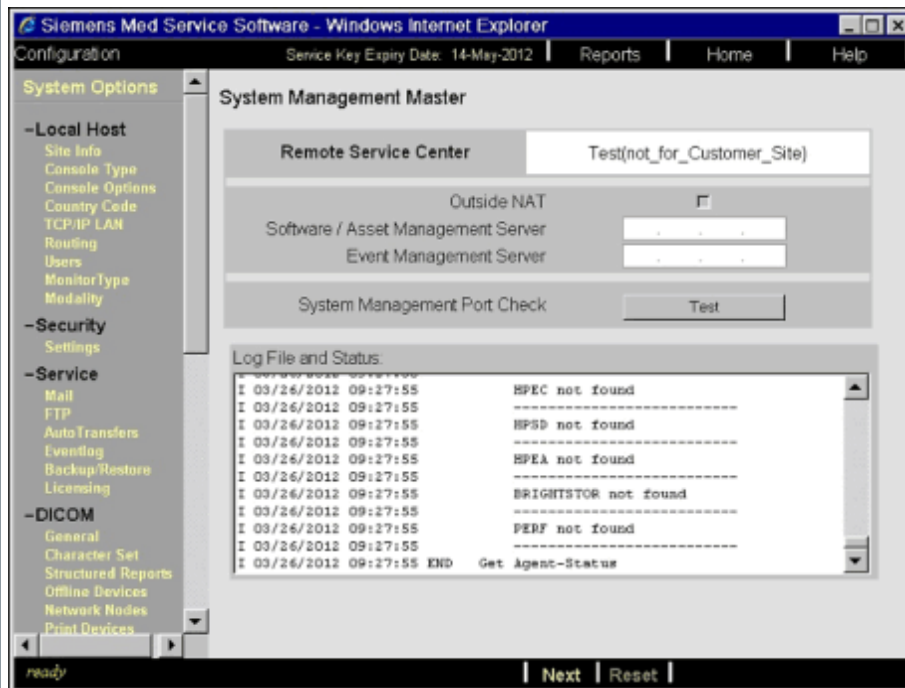

See the report file

C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html
for the Outside NAT configuration.

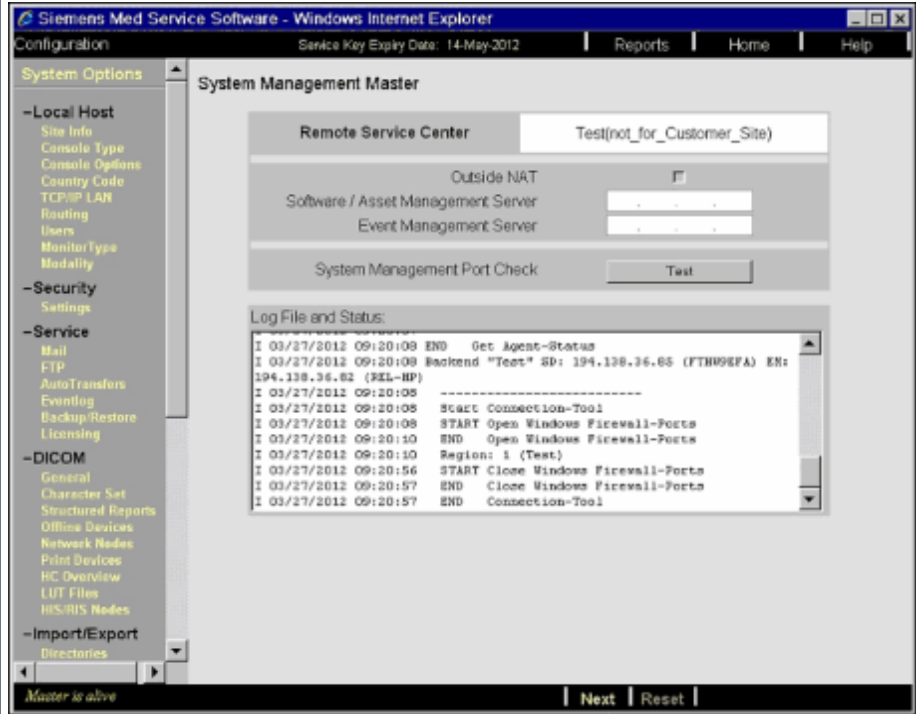ACTION	MRWP
■ Use the Test button to check whether the System Management Server can be accessed. Fig. 125: System Management master can be reached  The log file for the test can be found under the directory: C:\MedCom\service\mwtools\CCT\Logs	☐
In the event of problems with the connection test, contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	
■ Select Next.	☐


ACTION

MRWP

- Under **System Mgmt. Agent Controls**, set the following configuration:
Mark the checkboxes "Go Active" and "Go Online".
(The other checkboxes have already been marked).

☐

Fig. 126: System Management activation

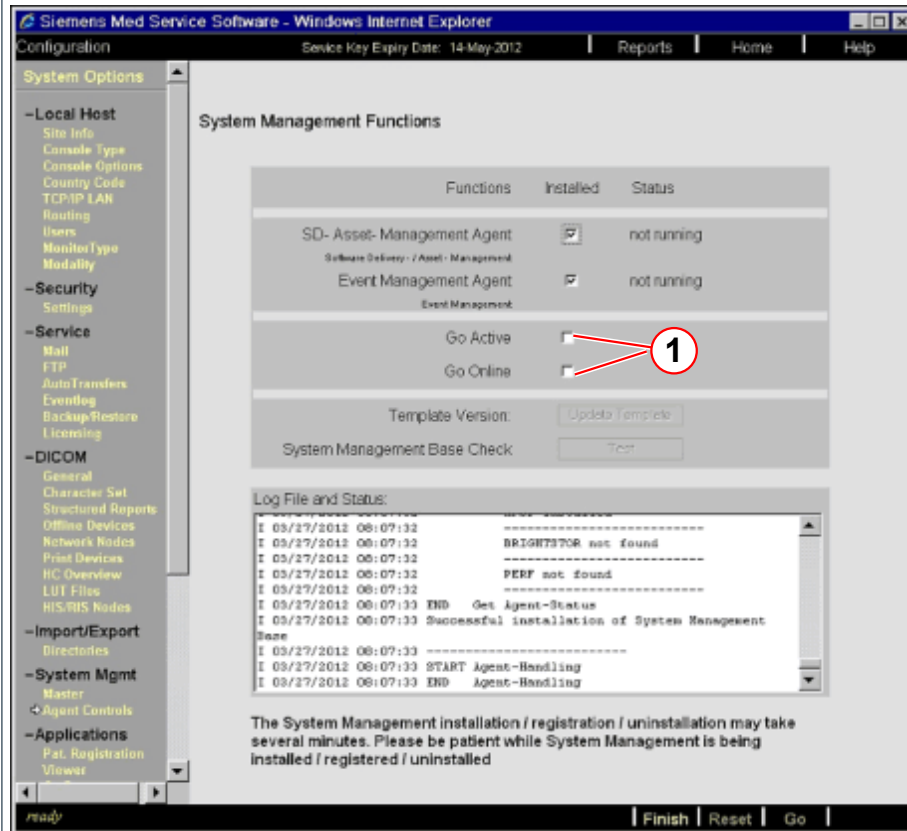



(1) Mark both checkboxes

In the Log File and Status area, the current status of System Management activation/installation is displayed.

- Click **Go**.

☐

Note: The installation/registration will take approximately 10 to 15 minutes.

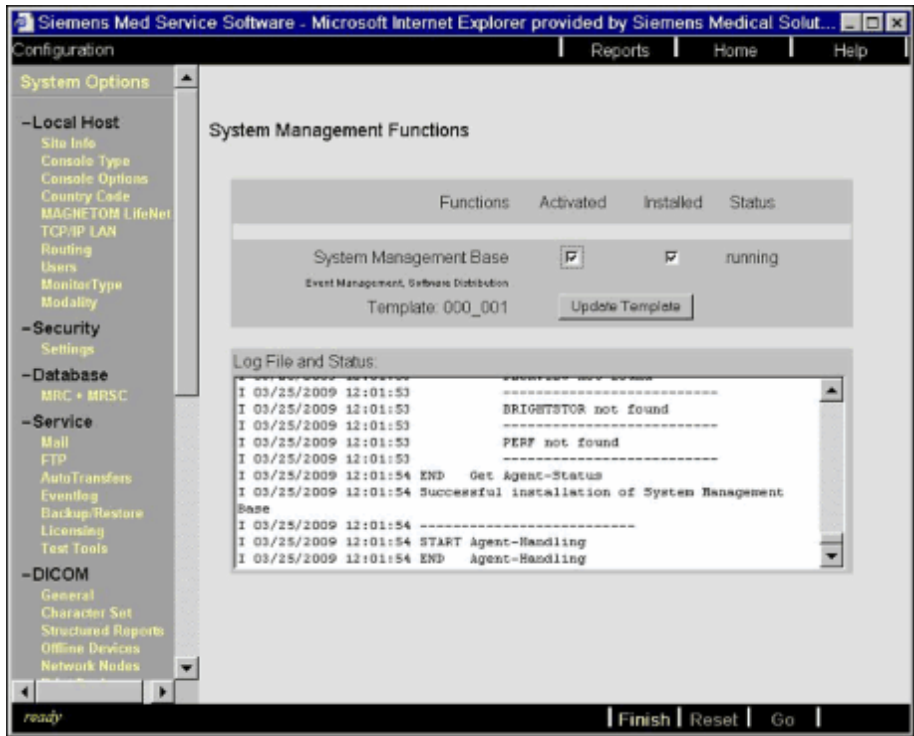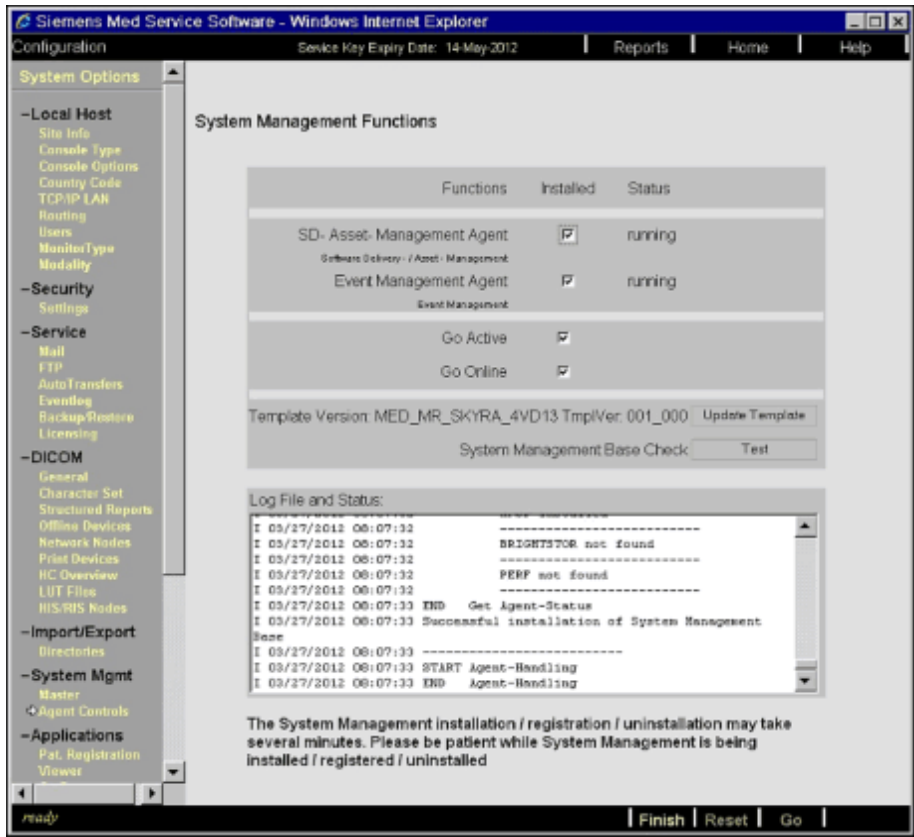ACTION	MRWP
■ If activation is successful, the following is displayed: Fig. 127: System Management function activation OK  Fig. 128: System Management activation OK  Select Finish to continue	☐

ACTION	MRWP
Click Finish and close the local service UI.	☐
In the event of problems (installation/activation fails a second time), contact the SRS help desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	

3.1.1.25 MRWP (MRSC) Enable the Auto Login Function



The following steps are required only if auto login is not enabled.

ACTION	MRWP
Select Options --> Service --> Local Service.	☐
Enter the valid service key and click OK.	☐
Select Configuration in the service platform.	☐
Select local host Users	☐
Mark the checkbox "Enable Auto login" to enable the auto login function.	☐
Click Save.	☐
Click Home to leave the configuration.	☐
Close the local service UI.	☐

3.1.1.26 MRWP (MRSC) Perform a COBR image backup on the MRWP (MRSC)



Please check the external media (e.g. USB disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.



A Win7 image backup with the script from the start menu is not recommended due to the limited size of the partition D.

ACTION	MRWP
- Create the COBR image backup of the system accordingly - Manual Backup of System Images to USB-Disk or - Manual Backup of System Images to DVD Note: Due to performance and time issues, it is recommended to use a USB disk for backup in place of DVD's.	☐

3.1.1.27 Deactivation of OLM (if applicable)



If the VE11A upgrade kit does already contain VE11A-SP01, the OLM has to be deactivated. If the VE11A upgrade kit does NOT contain VE11A-SP01, just VE11A, it needs to be checked if the OLM has already been deactivated either by UI MR019/15/R (VD13A-SP07) or by MR021/15/P (VD13C-SP03). If one of both has already been performed, it's mandatory that the VE11A-SP01 has to be installed afterwards (UI MR005/15/P) immediately.

- If the VE11A upgrade kit does contain VE11A only and not VE11A-SP01 (COBR label date lower than 24/07/2015) and the OLM was deactivated, it's mandatory that the VE11A-SP01 has to be installed afterwards (UI MR005/15/P).
If the OLM is still active, nothing else needs to be considered and you can proceed with "Final Check".
- If the VE11A upgrade kit does contain already VE11A-SP01 (COBR label date 24/07/2015 or above) and the OLM was already deactivated, nothing else needs to be considered and you can proceed with "Final Check".
If not (OLM still active) please follow the instructions according to Appendix (→ Deactivation of OLM / Page 216).

3.1.1.28 Final Check - to be performed for MRC and MRSC (if applicable)



If the customer uses a Network Time Protocol Server for synchronizing the time on the MR system, the configuration has to be done again. This configuration can not be done in the factory due to the check if the NTP server is reachable. (local service UI, Site Info, step 2)

- Verify the correct installation as follows:
In the main syngo MR UI, select
Help -> Info
Check whether the following entry is visible under the Siemens software package
In case it is a VE11A COBR-Recovery-Image DVD and has a creation date before 24/07/2015:
- N4_VE11A_LATEST_20140830_P3
In case it is a VE11A COBR-Recovery-Image DVD and has a creation date from 24/07/2015 onwards:
- N4_VE11A_LATEST_20140830_P12 (SP01)
or in case it is a VE11B COBR-Recovery-Image DVD
- N4_VE11B_LATEST_20150530_P1
The label
- N4_VE11A_LATEST_20140830_P3 / N4_VE11A_LATEST_20140830_P12 (SP01)
or
- N4_VE11B_LATEST_20150530_P1

are indicating a correct N4VE11A or N4VE11B base installation.

- Check the RDIAG configurations on the **MRC** and, if applicable, on the **MRSC** as well (the MRSC must also be connected to RDIAG) as follows:
 - Check the host name and IP address of the configured RDIAG server (Local Service -> Configuration -> Service Mail)

	Host name	IP address
Europe, Africa	LUX09505	194.138.39.18
North/South America	SRSACC1	129.73.116.92
Asia	sgpt806x	194.138.243.178

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- In the second step of the mail configuration, check the defined mail target (for the STE, the e-mail address must be changed according to the following table)

	Logical name for FTP target	Select the mail server host	E-mail address
Europe, Africa	SRS_mail	LUX09505	MRS MC.healthcare@siemens.com
	SaveLog	LUX09505	Savelog.healthcare@siemens.com
North/South America	SRS_mail	SRSACC1	MRUPTIME.healthcare@siemens.com
	SaveLog	SRSACC1	Savelog.healthcare@siemens.com
Asia	SRS_mail	sgpt806x	MRS MC.healthcare@siemens.com
	SaveLog	sgpt806x	Savelog.healthcare@siemens.com

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- In the second step of the FTP configuration (Local Service -> Configuration -> Service FTP), check the defined FTP target

	Logical name for FTP target	FTP host selection
Europe, Africa	SRS_ftp	LUX09505
North/South America	SRS_ftp	SRSACC1
Asia	SRS_ftp	sgpt806x

FTP user name: `rdiaftp`

Password: `siemens`

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Test the FTP connection to the RDIAG server as follows:

Local Service -> File&Image Tools

Source: File

Action: Transfer to System

File Type: Any file

From source: SRS_ftp

If a list from the directory /outgoing/mr (directory on the RDIAG server) is displayed, the connection and FTP access are OK

If the test is not OK, check the SRS configuration; refer to

- Click Home
- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **EL_Sample**. The following settings must be applied:

Active checkbox: must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **EL_Transfer**. The following settings must be applied:

Active checkbox: must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **IB_Sample**. The following settings must be applied:

The "Activated" checkbox must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **IB_Transfer**. The following settings must be applied:

The "Activated" checkbox must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **edb**. The following settings must be applied:

Active checkbox: must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **qdb**. The following settings must be applied:

Active checkbox: must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- The **UTILDATA** auto report is used for transferring the utilization data of the system. When the customer has a contract for utilization management, this auto report has to be marked as active also.



Prevent failure messages where the regional language setting did not match, e.g. UI language

- to guidance text within the workflow (Siemens pre-defined texts)
- to language text in error messages
- table position of the protocol after opening in a total different range and not useful for planning. (Due to mismatch between "." and ",")

See SKB0022441 entry and (→ Failure messages where the regional language setting did not match, e.g. UI language / Page 182).

3.1.1.29 Customer Information

The mandatory (5 days on-site) Appl. Training is being provided via (→ Application Training Order Guide,

<https://intranet.healthcare.siemens.com/cms/CS/en/application/MR/Pages/Default.aspx>).

CS-Homepage - For Application - MR - Training Order Guide MR

Potentially additional training has to be determined individually.

3.1.1.30 Final Work Steps

1. For systems equipped with the Expert-i license:

Follow the instructions in document MR00-000.814.28 to update existing client-PCs using the Expert-i client software delivered with the CSSI upgrade.

Otherwise scrap the Expert-i client software.

3.2 CSSI Upgrade Second Steps for VE11C




After receiving the CSSI kit for the software upgrade to VE11C, the following steps have to be performed.

3.2.1 Work Steps

3.2.1.1 Overview

Perform the CSSI upgrade

Fig. 129: CSSI Final Steps Flowchart

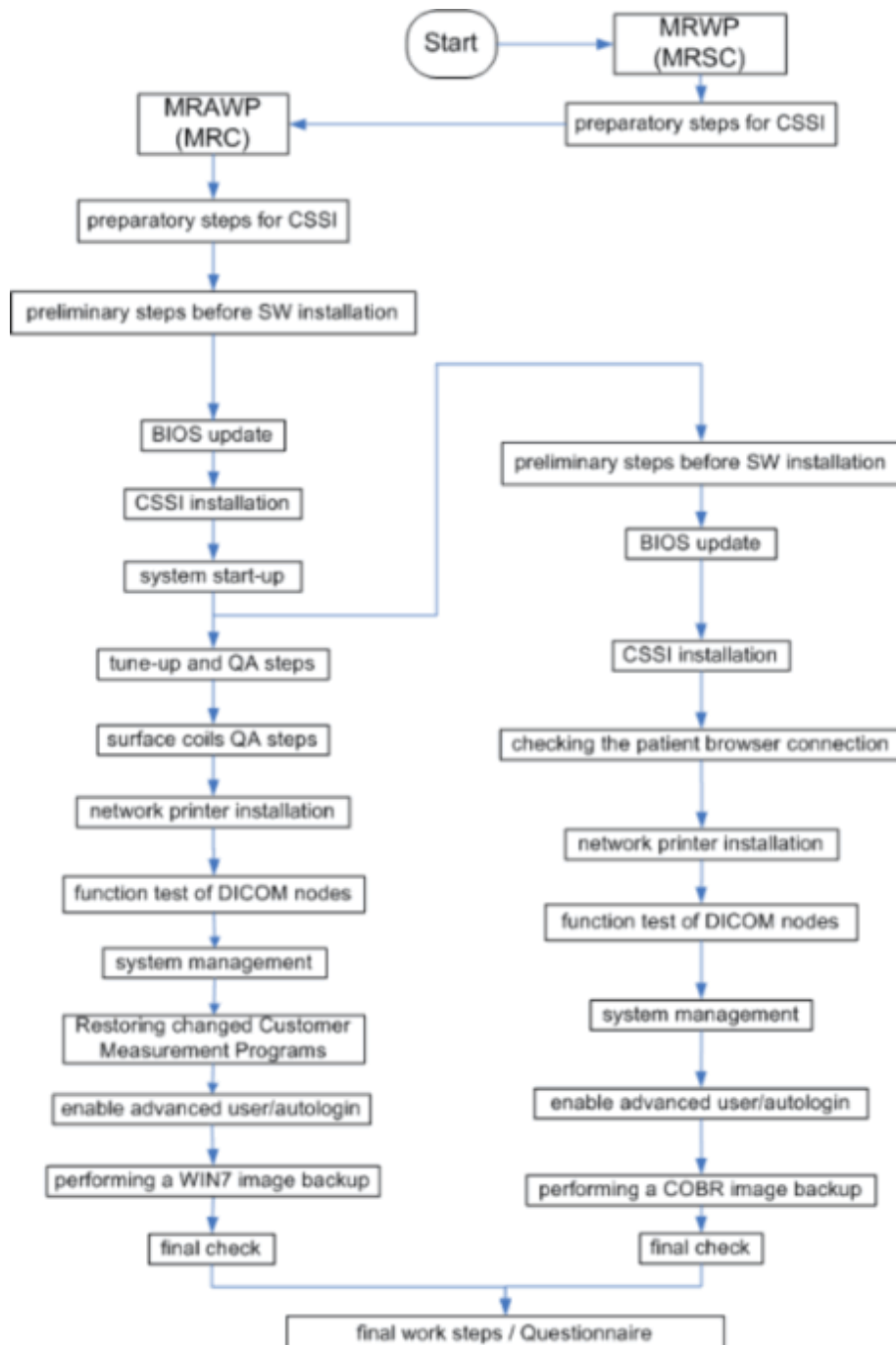

3.2.1.2 MRWP (MRSC) Preparatory Steps for CSSI Upgrade with the Old Software Version



If you have a satellite console with a shared database connected to the main console, perform the following steps first.

If no MRWP is present, skip the following steps and start at (→ MRAWP (MRC) Preparatory Steps for CSSI Update with the Old Software Version / Page 118)



Please check if any local or network printer is configured on MRWP and clarify with the customer representative the availability of installation material (Software, Installation Manual, etc.) if necessary. Please note the configuration data.

Performing a Master Backup on the MRWP (MRSC)



Please check if any local or network printer is configured on MRC and clarify with the customer representative the availability of installation material (Software, Installation Manual, etc.) if necessary. Please note the configuration data.



If required (configuration is changed) the step (→ Collecting MRWP Configuration Data from Local Service / Page 33) can be done too.

ACTION	MRWP
Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Master Backup VD13 or VE11A or VE11B MRWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK.	☐

ACTION	MRWP
Note: ■ Check the output log for errors. If any other error occurs apart from those mentioned below, the master backup has to be created again. ■ In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed" appears. Click "OK" to continue ■ If the following error messages are displayed at the end of the backup routine, ignore them. The backup is completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
■ Remove the CD.	☐

MNP Deregistration on the MRWP



In order to avoid problems with the MNP registration with the new Numaris version, it is recommended to perform the deregistration with the old Numaris version.

ACTION	MRWP
■ In the local service UI, select Configuration.	☐
■ Under System Mgmt. Agent Controls, set the following configuration: Remove the mark from the checkboxes "Go Active" and "Go Online".	☐
■ Click Go.	☐
The deregistration of the MNP is done. This will take about 5 minutes.	
■ Click Finish.	☐
■ Click Home and close the local service UI.	☐
■ Shut down the MRWP.	☐
■ Switch off the MRWP.	☐

3.2.1.3 MRAWP (MRC) Preparatory Steps for CSSI Update with the Old Software Version

- 

If you have a satellite console with a shared database connected to the main console, start with the action item:
(→ MRWP (MRSC) Preparatory Steps for CSSI Upgrade with the Old Software Version / Page 116)
This has to be done first. If this step has already been completed on the MRWP, continue with the MRAWP.
- 

Please check if any local or network printer is configured on MRAWP and clarify with the customer representative the availability of installation material (Software, Installation Manual, etc.) if necessary. Please note the configuration data.

Performing a Master Backup on the MRAWP (MRC), for changed Customer Protocols

- 

Please check if any local or network printer is configured on MRC and clarify with the customer representative the availability of installation material (Software, Installation Manual, etc.) if necessary. Please note the configuration data.
- 

If required (configuration is changed) the step (→ Collecting MRAWP Configuration Data from Local Service / Page 35) can be done too.

ACTION	MRAWP
■ Perform a Master Backup in local service. The Master Backup includes the Numaris4 Package, Software Settings, and Security Settings. Label an empty CD-R medium with Master Backup VD13 or VE11A or VE11B MRAWP (actual date) Start local service, select Backup&Restore, insert the labeled CD-R medium into the [CD-R] drive (S:\), select Backup in the Command menu, select drive S:[CD-R] in the Drives menu, select Master-Backup in the Package menu, start the backup by clicking OK.	☐

ACTION	MRAWP
Note: ■ Check the output log for errors. If any other error occurs apart from those mentioned below, the master backup has to be created again. ■ In some cases, the following popup message: "Please insert the media into Drive S:\ and press 'OK' to proceed" appears. Click "OK" to continue ■ If the following error messages are displayed at the end of the backup routine, ignore them. The backup is completed correctly. - Error Execute starting components... 3 rti.exe runComp an_server rti->runComp(an_server, Pid, Container) failed and/or - Error .. Execute starting components... 2 rti.exe runComp stu_server rti->runComp(stu_server, Pid, Container) failed .. and/or - Error Execute starting components... 1 rti.exe runComp ima_server rti->runComp(ima_server, Pid, Container) failed ...	
■ Remove the CD.	☐

MNP Deregistration on the MRAWP



In order to avoid problems with the MNP registration with the new Numaris version, it is recommended to perform the deregistration with the old Numaris version.

ACTION	MRAWP
■ In the local service UI, select Configuration.	☐
■ Under System Mgmt. Agent Controls, set the following configuration: Remove the mark from the checkboxes "Go Active" and "Go Online".	☐
■ Click Go.	☐
The deregistration of the MNP is done. This will take about 5 minutes.	
■ Click Finish.	☐
■ Click Home and close the local service UI.	☐

Manual backup of voice commands on MRAWP



Due to a bug in the backup program, it is recommended to copy the voice command files to an external medium and then to restore these in the new software version.

Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

ACTION	MRAWP
■ Enable the advanced user - Run the Advanced User program from the Start menu (press <CTRL> + <ESC> to access the Start menu) - From the Start menu, select All Programs -> Advanced User. - In the popup displayed, enter the advanced user password and press the OK button - Confirm Advanced User enabled popup with OK.	☐
■ Connect a USB medium.	☐
■ Start Windows Explorer: - Press <CTRL> + <ESC> to access the Start menu - In the Start menu select All Programs ⇒ Accessories ⇒ Windows Explorer Windows Explorer starts.	☐
■ Navigate to the directory C:\MedCom\MriSiteData	☐
■ Copy the subdirectory PIA and its contents to the USB medium. The subdirectory PIA contains the voice command files.	☐
■ Close Windows Explorer.	☐
■ Remove the USB medium.	☐

3.2.1.4 MRAWP (MRC) Preliminary Steps before the Software Installation

ACTION	MRAWP
■ Ensure that the network cables at the back of the console are connected. - MRC: Two network cables at the back of the computer.	☐
■ Ensure that only the monitor and the PDD are connected. Disconnect any other connected USB customer devices (such as USB hard disks). Keep only the USB mouse, the USB keyboard and powerloading stations for ECG connected.	☐
■ If applicable, check whether the connected local printer (connected to the parallel port) is switched ON. During the installation, the printer is automatically detected and the correct printer driver is installed.	☐
■ Disable the Autologin function in the local service UI.	☐

ACTION	MRAWP
■ Perform a shutdown and switch off the system on the alarm box.	☐
■ Remove the power cord connector from the MARS in the Electronic Power Cabinet. (This has to be done because otherwise the MARS installation can run into problems)	☐



If you have the following old Numaris version VD13* (except VD13D) and a HPz420 host, the CSSI VE11C update kit additionally includes a memory upgrade to 32 GBytes. Please perform the step (→ Memory upgrade to 32 GByte / Page 241)


The Prisma and PrismaFIT systems with VD13D and HPz420 already have the 32 GByte memory installed.


3.2.1.5 MRAWP (MRC) BIOS Update









Ensure that the power cord connector is removed from the MARS.




ACTION	MRAWP
■ Switch on the system on the alarm box	☐
■ Login as user administrator	☐
■ Insert the first DVD of the supplied CSS image set labeled COBR Recovery Image MR system software VE11C MR system serial no: <x> AWP into the CD/DVD burner drive (S:\).	☐
■ Insert the supplied CD labeled with "MR-VD BIOS CD VE11.." into the host CD/DVD-ROM Drive (R:\)	☐
■ Invoke a command shell (Start ⇒ All Programs ⇒ Accessories ⇒ Command Prompt)	☐
■ In the command shell, type in: S:\Start.bat	☐

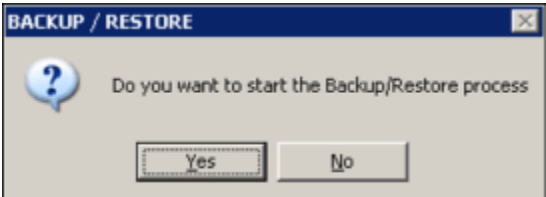
ACTION	MRAWP
Confirm the Installation popup with "Yes" <i>Fig. 130: Backup/Restore popup</i> 	☐
The running script performs a reboot at the end. This may take some time. Please wait until this is finished.	☐
The BIOS CD will be booted automatically.	☐
The BIOS update runs automatically (no user action is necessary) and ends when the following message is displayed: <pre>BIOS Installation Done - System Ready for Software Installation Please remove BIOS CD ! Please Turn Off Power or Reset System ! C:\> The BIOS update takes about 5 minutes.</pre>	☐
Remove the supplied BIOS CD from the CD/DVD drive (R:\)	☐
Remove the first DVD of the supplied CSS image set labeled COBR Recovery Image MR system software VE11C MR system serial no: <x> AWP from the CD/DVD burner drive (S:\).	☐
Insert the first DVD of the supplied CSSI image set labeled COBR Recovery Image MR system software VE11C MR system serial no: <x> AWP into the CD/DVD-ROM drive (R:\).	☐
Press key combination <CTRL> + <ALT> + (on a German keyboard use <Strg> + <Alt> + <Entf>) to reboot and continue with "CSSI installation MRC".	☐

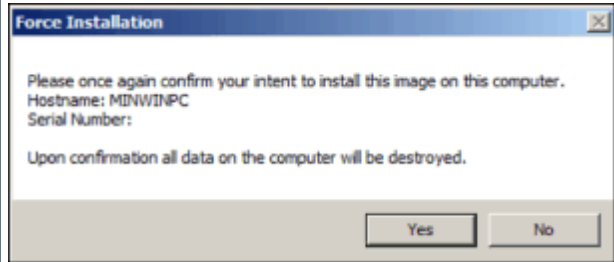
ACTION	MRAWP
Confirm the Force Installation popup (example is shown below) with Yes. <i>Fig. 133: Force Installation popup</i> A Windows-style dialog box titled "Force Installation". The text inside says: "Please once again confirm your intent to install this image on this computer. Hostname: MINWINPC Serial Number: Upon confirmation all data on the computer will be destroyed." At the bottom are "Yes" and "No" buttons. This Force Installation popup is displayed twice. Confirm it each time with Yes.	☐
The system restore process starts now. The following shows an example of the screen display. <i>Fig. 134: COBR restore process start</i> The image shows a computer screen during a system restore process. The background is a blue Siemens Healthcare Sector logo with the text "COBR VA31C". In the foreground, a black terminal window displays white text showing the progress of the restore process, including file copying and partition creation.	☐

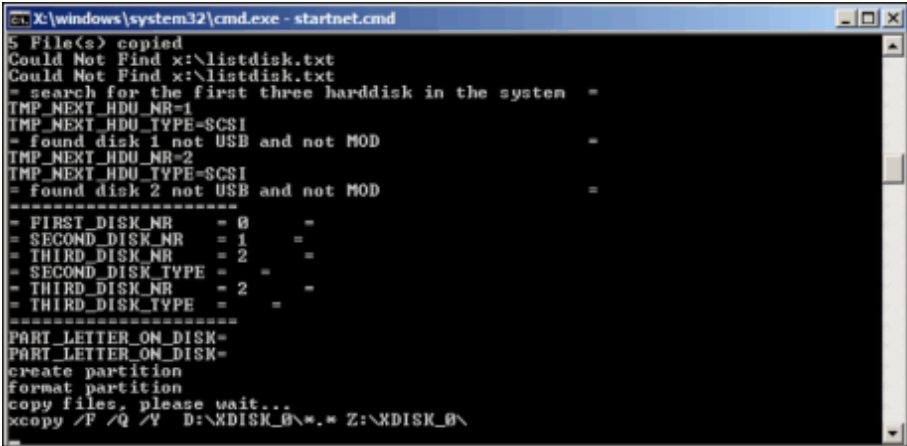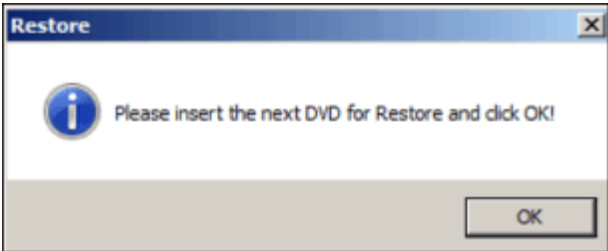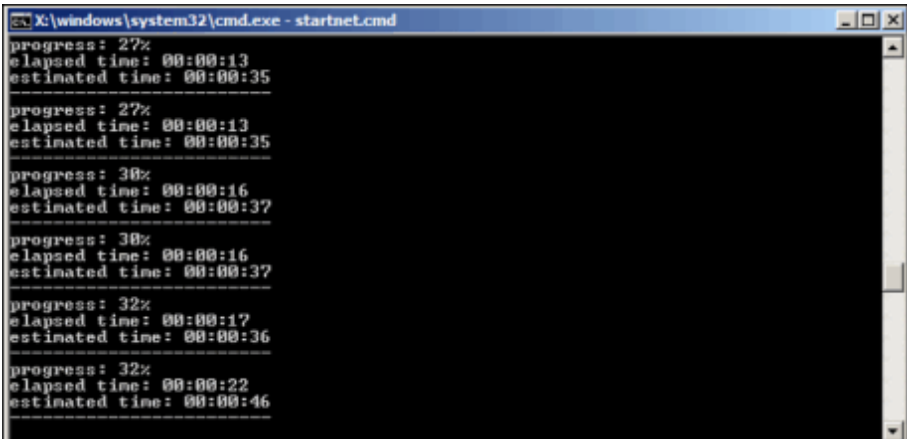
ACTION	MRAWP
First all image files from the DVD for the system disk (disk1) are copied from the DVD to the hard disk. Fig. 135: Copy image files  After the files are copied, the DVD is automatically ejected.	☐
Remove the ejected medium from the CD/DVD-ROM drive, insert the next DVD from the system image DVD set, and close the drive.	☐
Confirm the request popup with OK. Fig. 136: Next DVD request popup 	☐
The copy procedure continues. The next DVD will be requested automatically until all required data has been copied. The copy and the confirmation of the popup (→ Fig. 136 Page 125) have to be repeated for all available CSSI DVDs.	
After the last CSSI DVD has been copied, the Image Restore starts Fig. 137: Image restore progress output 	☐

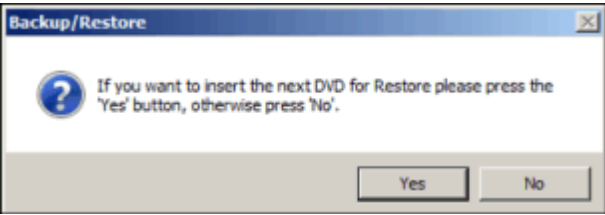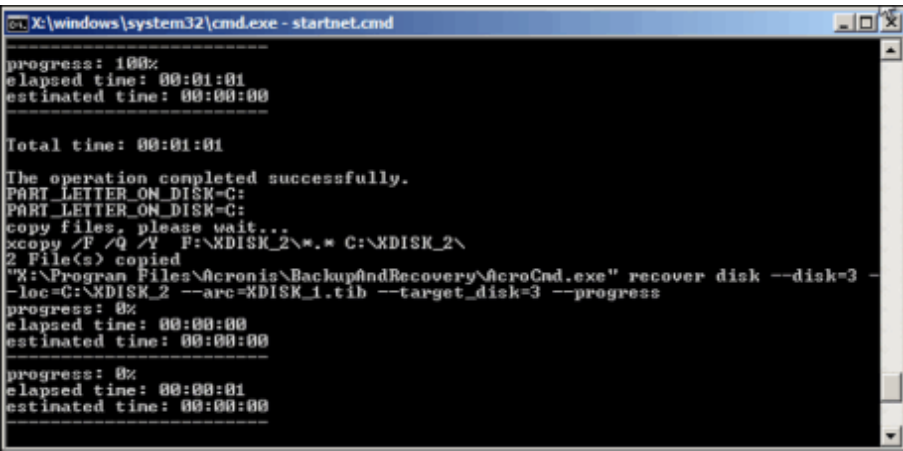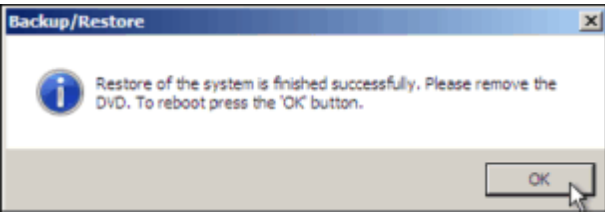ACTION	MRAWP
Wait until the following popup appears. <i>Fig. 138: Next DVD request popup</i> 	☐
Confirm the "Next DVD request popup" shown above with "No".	☐
Now the last disk is restored from the image. <i>Fig. 139: Last disk restore</i> 	☐
When all data has been successfully restored, the following popup is displayed. <i>Fig. 140: Backup/Restore successfully finished popup</i> 	☐
Remove the ejected medium, close the drive, and confirm this popup with OK.	☐
The system is rebooted.	☐
Note: The following behavior is new in the VE11C CSSI update.	
After the reboot, the autologin fails with an error message. Please wait until a further reboot is done (takes about 30 seconds). The next autologin will be successful. Note: Please don't click or hit any key - please just wait.	☐

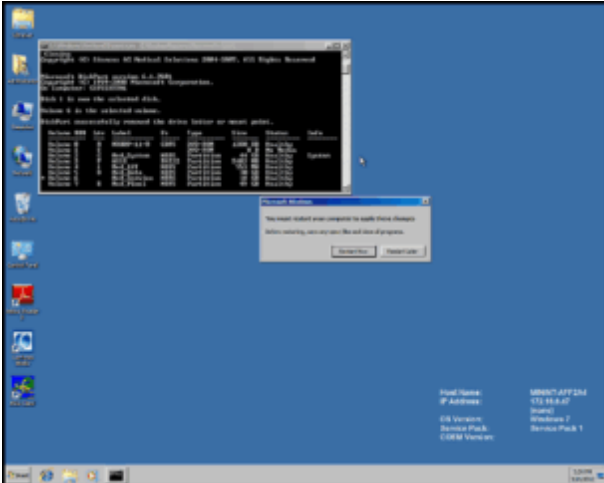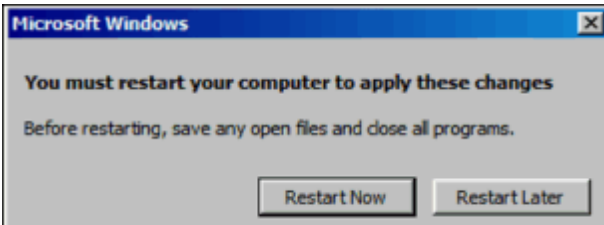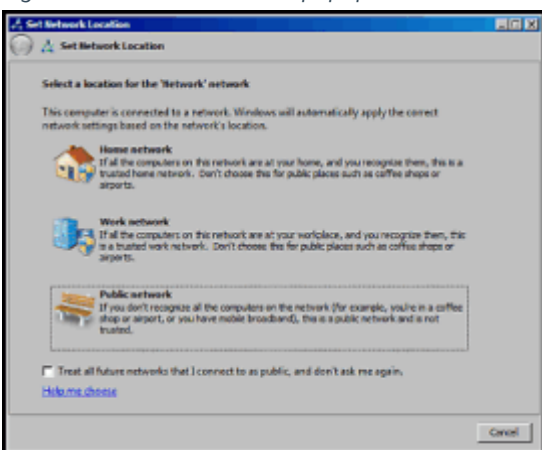ACTION	MRAWP
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 10 minutes. Note: Don't accept any displayed popup or press a button on any displayed popup.	☐
During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed: <i>Fig. 141: Displayed popup</i>  Note: Don't accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don't touch or accept it. Keep it displayed on the screen. It will be closed automatically. <i>Fig. 142: Windows popup</i>  Note: The following popup can be ignored too, if it is displayed. <i>Fig. 143: Network location popup</i> 	

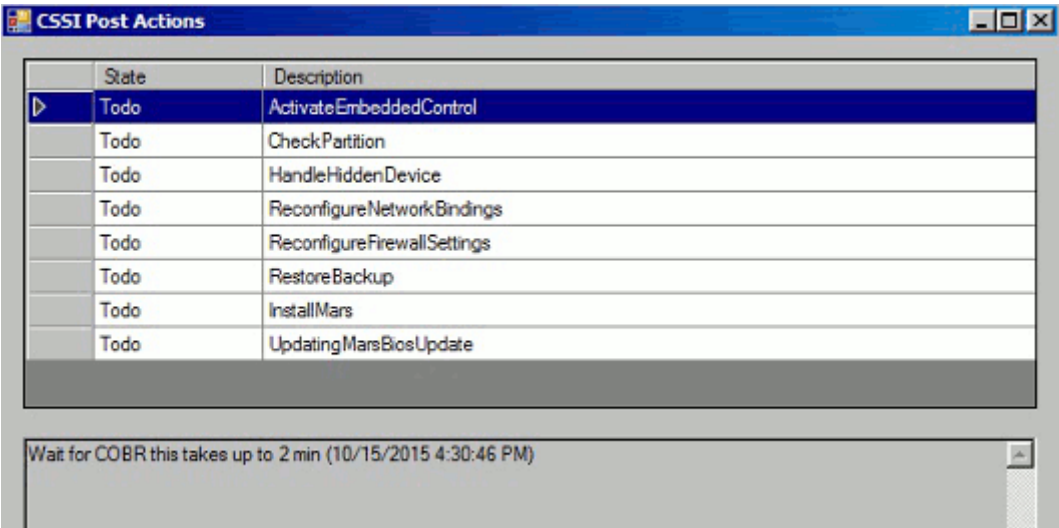ACTION		MRAWP
The CSSI Post popup is displayed (see)(→ Fig. 144 Page 128)		
Fig. 144: CSSI Post Action popup		
		
The following CSSI Post Action steps are performed:		
■ ActivateEmbeddedControl (duration under 1 min.) ■ CheckPartition (duration under 1 min.) ■ HandleHiddenDevice (duration under 1 min.) ■ Reconfigure network connections (duration under 1 min.)		
■ If hidden network devices are found and removed, the system is rebooted automatically. After the autologin as syadmin the CSSI Post Actions popup (→ Fig. 144 Page 128) is displayed again.		☐
■ The following CSSI Post Action steps are performed: - Reconfigure firewall settings (duration under 1 min.) - Restore backup (only IP settings, duration about 1 min.)		☐
■ The system is rebooted automatically.		☐

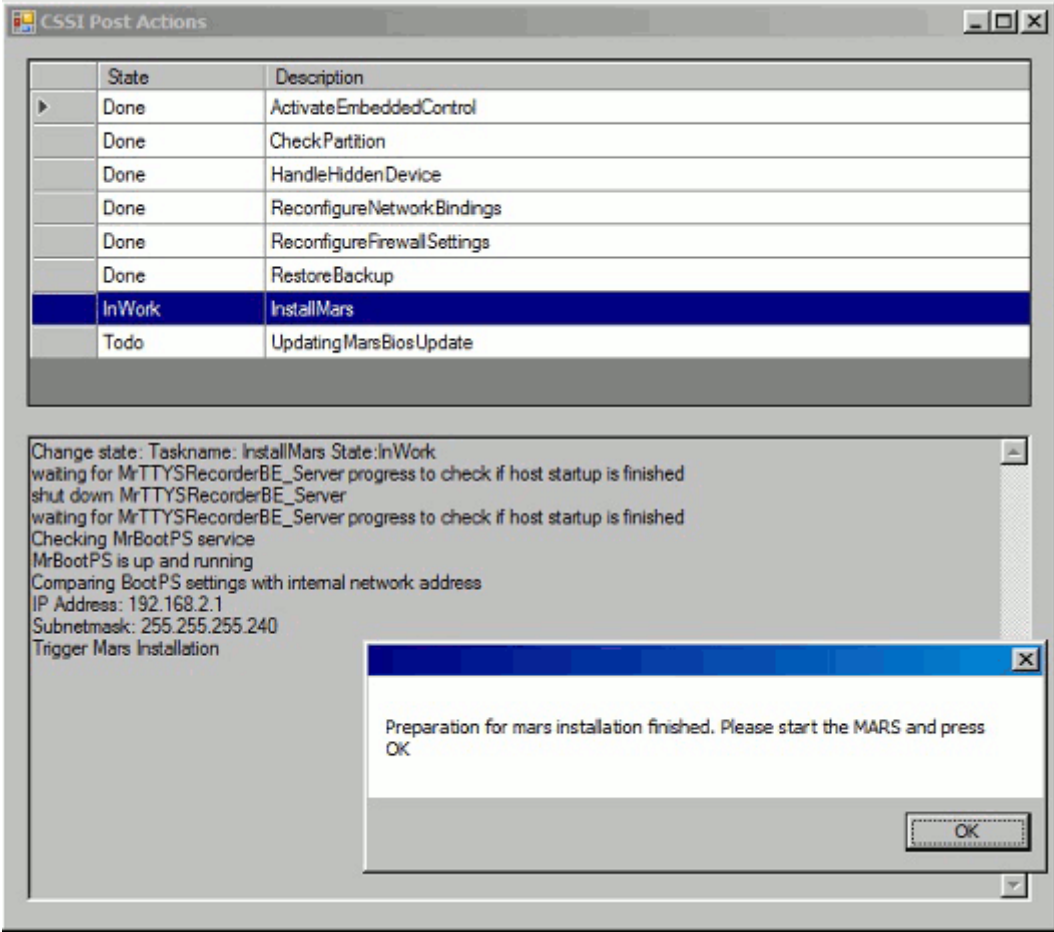ACTION	MRAWP
■ An autologin as user "syadmin" is performed and the CSSI Post Actions popup (→ Fig. 144 Page 128) is displayed. The following CSSI Post Action steps are performed: - Checking prerequisites for MARS BIOS update Wait until you get the Please start the MARS... popup (takes about 5 min.). <i>Fig. 145: MARS Start popup</i> 	☐
■ Connect the power cord connector from the MARS in the Electronic Power Cabinet.	☐
■ Confirm the Please start the MARS... popup with OK.	☐

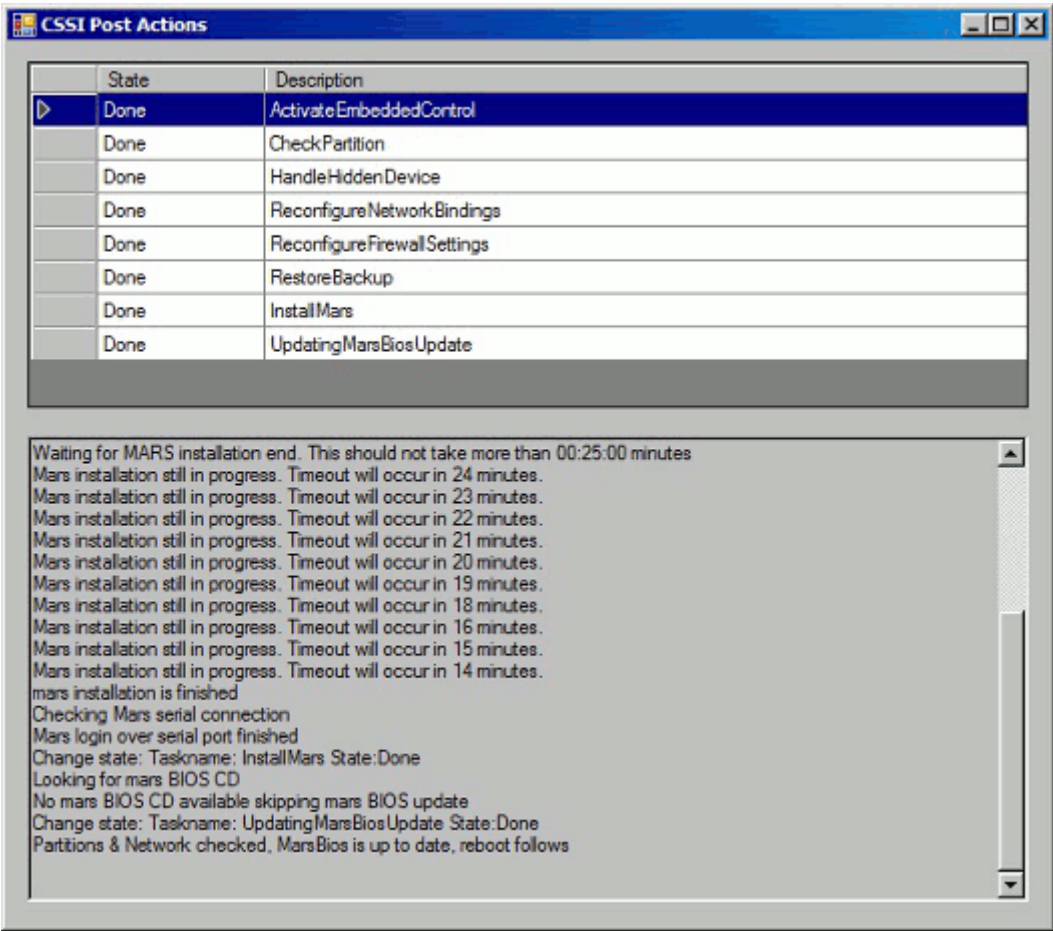ACTION	MRAWP
■ The MARS is booted up and the following CSSI post action step continues: - Install MARS (duration about 15 min.) A reboot is performed automatically after this step is done. Fig. 146: MARS install OK reboot follows 	☐

ACTION	MRAWP
Note: In case of an error, the CSSI Post Action stops at the step which shows its status as failed. The reboot is then not done automatically. In this case, please perform the following workaround. Run the CSSI Post Actions a second time. - Close the CSSI Post Actions by pressing the X. - Start a cmd tool window (Start ⇒ All Programs⇒ Command Prompt). - In the cmd tool window change the directory to <code>C:\MedCom\utils</code> - Execute the CSSI Post actions again by typing in the command: <code>CSSIPostRestoreProcess.exe ↵</code> Now the CSSI Post actions tool starts again. This takes about 20 minutes. After this time a reboot is done automatically. If the CSSI Post actions tool fails again, please contact the HSC MR in Erlangen.	
- Continue with one of the following steps, depending on the old Numaris version installed: - If old version was VE11A/VE11B, continue with (→ Worksteps for the update from VE11A/VE11B / Page 131) - If old version was VD13, continue with (→ Worksteps for the update from VD13 / Page 133)	☐

3.2.1.7 MRAWP (MRC) System Start-up




If the customer has a TIM docking table (mobile), the table must be connected during the whole firmware/loadware update process.

If the TIM docking table is removed during the firmware/loadware update process, the PSY (Power Stage Y) will be destroyed.


Worksteps for the update from VE11A/VE11B

ACTION	MRAWP
Wait until the main UI is displayed.	☐

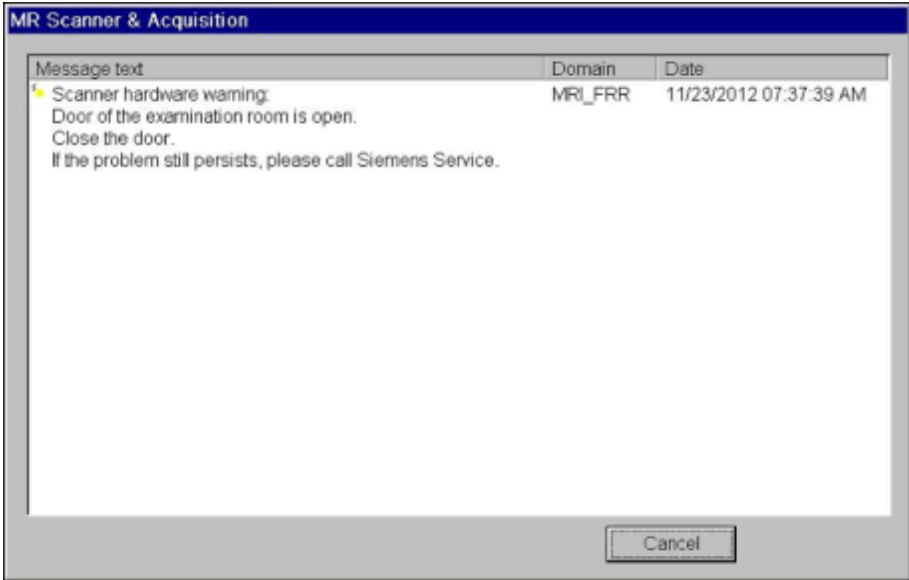ACTION	MRAWP
■ The firmware update starts, and the MR Scanner & Acquisition message popup is automatically displayed. Fig. 147: Example for MR Scanner & Acquisition popup 	☐

Note: A message is displayed for each component when the download is started. This message also includes the duration of the download. An example is shown below:

Fig. 148: Message popup

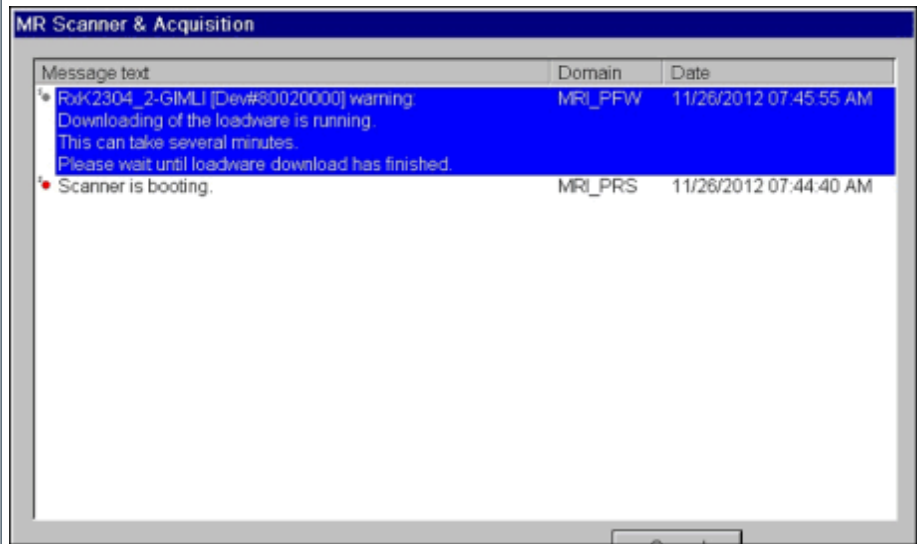

This message will disappear automatically when the component download is finished. -> You can check if an update has been running recently.

Note: For the PTAB, an error will be displayed until the firmware has been loaded successfully.


If no message is displayed, the wrong language setting has been used. Please correct it accordingly (→ Correct syngoUI Language Settings / MR00-000.814.55) after the firmware update. In this case please observe the firmware update in the Windows event log.

ACTION	MRAWP
Wait until the component firmware update has finished and the scanner is ready for measurement. This can take up to 30 minutes max.	☐
Continue with (→ Worksteps for the update from VE11A/VE11B / Page 136)	☐



During the firmware update of the SEP chiller, the refrigerator (i.e. helium compressor) is switched off. It will be automatically switched on again after about 20 minutes.

Worksteps for the update from VD13

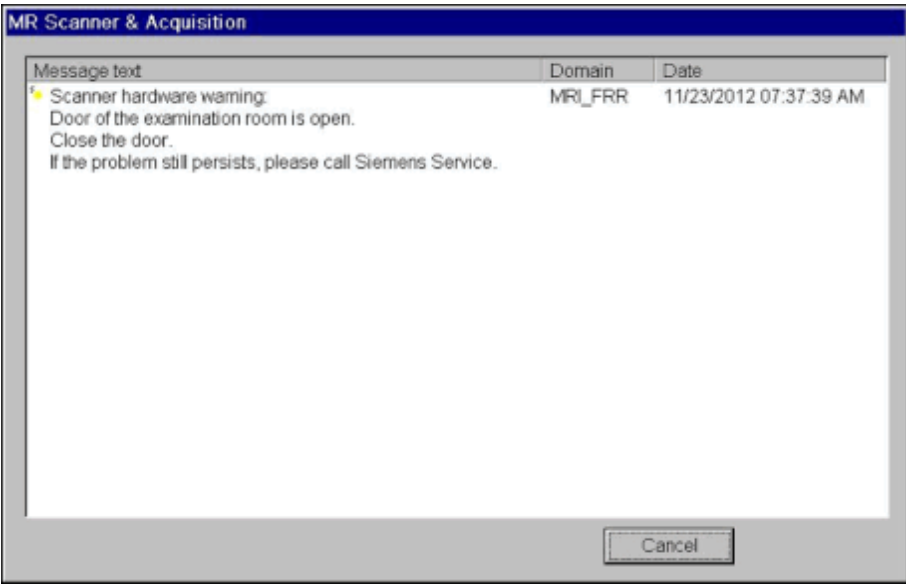
ACTION	MRAWP
Wait until the main UI is displayed.	☐
The firmware update starts, and the MR Scanner & Acquisition message popup is automatically displayed. Fig. 149: Example for MR Scanner & Acquisition popup 	☐

ACTION

MRAWP

Note: A message is displayed for each component when the download is started. This message also includes the duration of the download. An example is shown below:

Fig. 150: Message popup

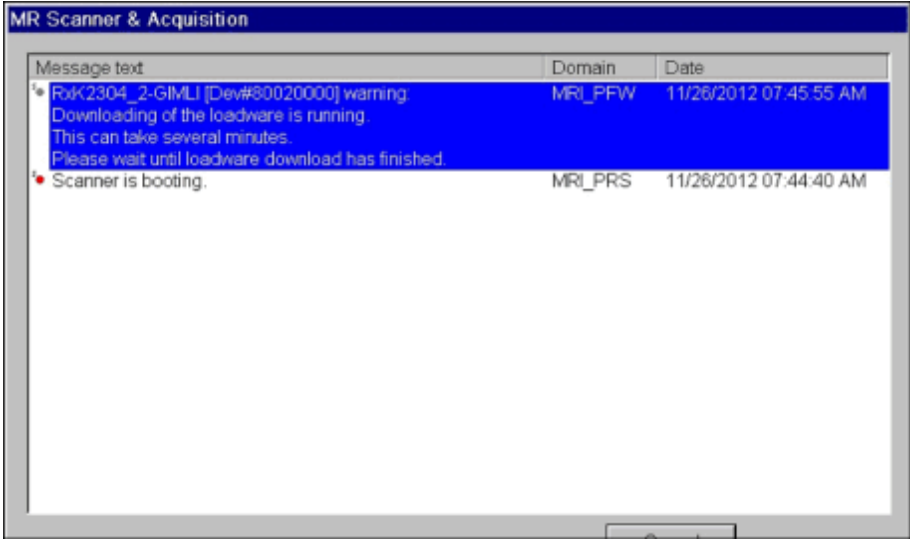


This message will disappear automatically when the component download is finished. -> You can check if an update has been running recently.

Note: For the PTAB, an error will be displayed until the firmware has been loaded successfully.

If no message is displayed, the wrong language setting has been used. Please correct it accordingly (→ Correct syngoUI Language Settings / MR00-000.814.55) after the firmware update. In this case please observe the firmware update in the Windows event log.

The firmware update time duration is about 60 minutes.

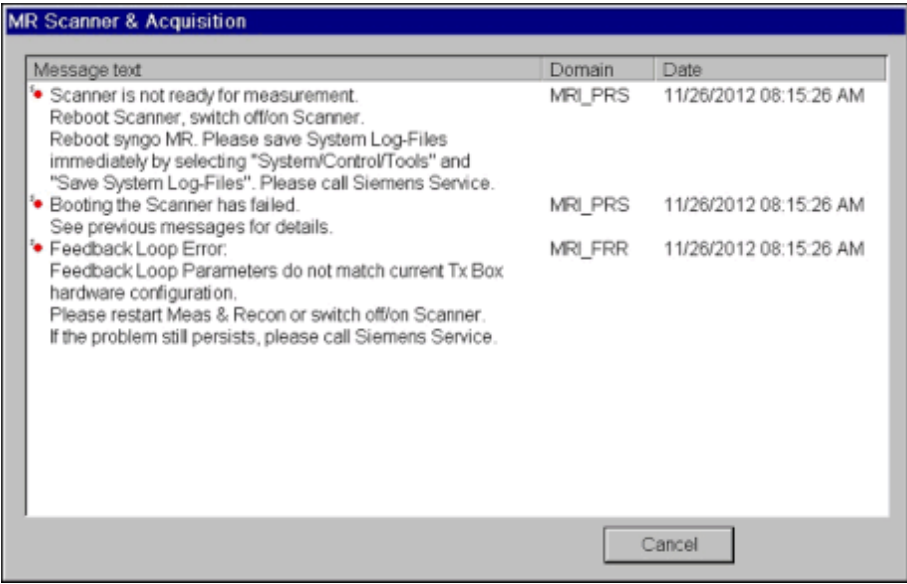
ACTION	MRAWP
■ Wait until the component firmware update has finished. The firmware update has finished when the following message is displayed in the MR Scanner & Acquisition message popup: - Feedback Loop Error: Feedback Loop Parameters do not match the current TX-Box hardware configuration. Please restart Meas & Recon or switch off/on the Scanner. If the problem persists, please call Siemens Service. The reason for this error message is the new Tx-Box firmware which resets the parameters. It can take up to 60 minutes max. until this message appears. <i>Fig. 151: Message popup "Feedback Loop Error"</i> 	☐
■ Continue with (→ Worksteps for the update from VD13 / Page 137)	☐

3.2.1.8 MRAWP (MRC) Tune-up and QA Steps

Worksteps for the update from VE11A/VE11B

ACTION	MRAWP
Restore the Tuneup status: Go into local service UI: - Select Utilities ⇒ Escape to OS and type in the command <code>MrSesoRestoreTuneUpStatus.bat</code> and press <ENTER> Fig. 152: SESO output  This step restores the last TuneUp&QA steps. The QA steps for the surface coils are restored too.	☐
Click Home and select Tune Up.	☐
Perform all tune-up steps that have the status "ToDo".	☐
Note: If you have a satellite console MRWP (MRSC), you can now start with the step (→ MRWP (MRSC) Startup / MR00-000.814.55) in parallel with the tune up and QA.	
- Perform all QA steps that have the status "ToDo". If all QA steps are set to "Done" perform the following steps: - BC Power Losses Check - BC Image Check - Image Orientation Check	☐
Continue with (→ MRAWP (MRC) Surface Coils QA Steps / Page 138)	☐

Worksteps for the update from VD13

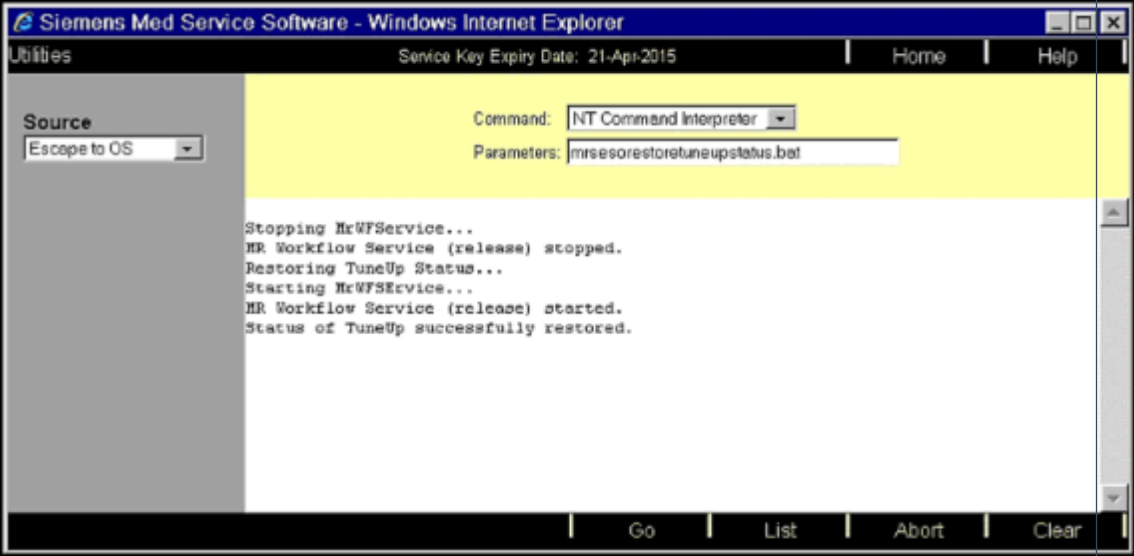ACTION	MRAWP
■ Restore the Tuneup status: Go into local service UI: - Select Utilities ⇒ Escape to OS and type in the command <code>MrSesoRestoreTuneUpStatus.bat</code> and press <ENTER> <i>Fig. 153: SESO output</i>  This step restores the last TuneUp&QA steps. The QA steps for the surface coils are restored too.	☐
■ Click Home and select Tune Up.	☐
■ Perform the following tune-up step first (keep the default selection of the tune-up steps): - Feedback Loop Calibration During this step the MARS is rebooted.	☐
■ Perform all other tune-up steps that have the status "ToDo"	☐
Note: If you have a satellite console MRWP (MRSC), you can now start with the step (→ MRWP (MRSC) Startup / MR00-000.814.55) in parallel with the tune up and QA.	
■ Perform all QA steps that have the status "ToDo"	☐
■ Continue with (→ MRAWP (MRC) Surface Coils QA Steps / Page 138)	☐

3.2.1.9 MRAWP (MRC) Surface Coils QA Steps



For all coils used before the software update, it is not necessary to perform a Coil QA. The QA status is restored.


3.2.1.10 MRAWP (MRC) Network Printer Installation


ACTION	MRAWP
Install the network paper printer according to the software installation manual (this step is not part of the dynamic data) if previously installed on the MRAWP or MRWP.	☐
If applicable, set the installed printer as the default printer. Select Start -> Settings -> Printers and Faxes Select the printer and click the right mouse button. In the menu displayed, select the "Default printer" entry.	☐

3.2.1.11 MRAWP (MRC) Configuration checks

ACTION	MRAWP
Select Options ⇒ Service ⇒ Local Service.	☐
Enter the valid service key and click OK.	☐
Select Configuration in the service platform.	☐

ACTION	MRAWP
■ Select Monitor Type under Local Host and check the settings as follows: - Monitor Type: Color and Siemens (True Color) - Multi-Monitor settings Only "Start syngo in single monitor" is marked Fig. 154: Monitor settings SESO 	☐
■ Click Save if changes have been made and Finish	☐

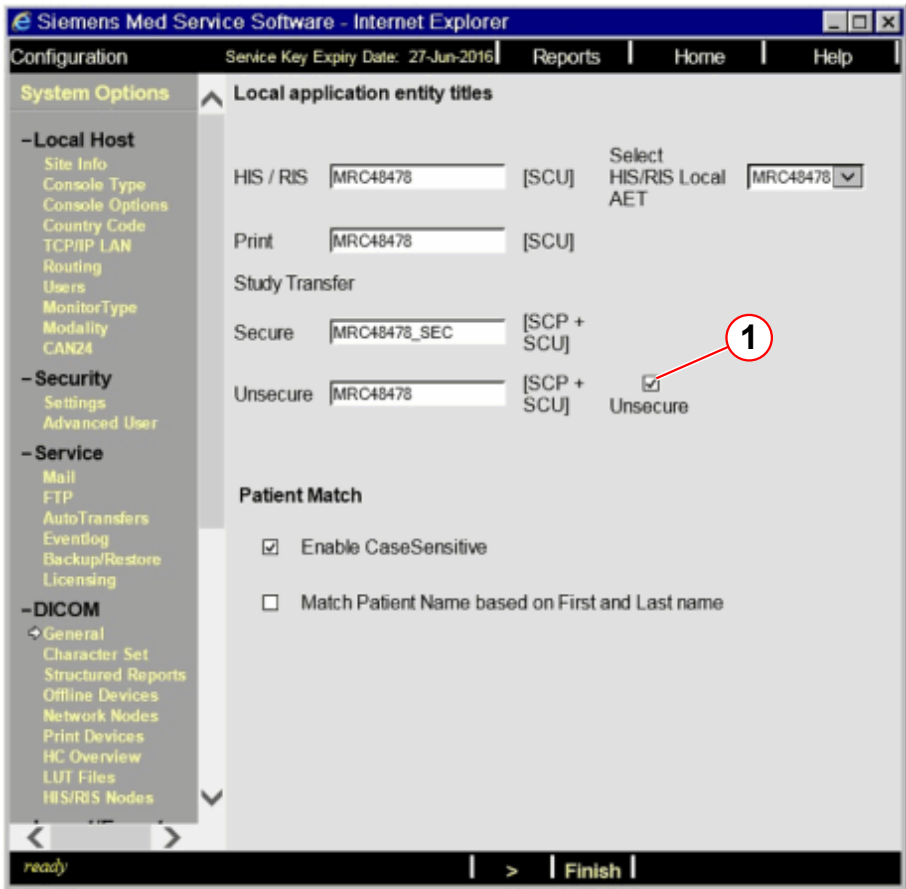ACTION		MRAWP
■ Under DICOM, select General and check if the checkbox Unsecure is marked. Fig. 155: Example form sim console 		☐
(1) Unsecure checkbox must be marked		
■ Click Save if changes have been made and Finish		☐
■ Under DICOM, select Network Nodes and press under Host Properties >.		☐
■ On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unsecure. Fig. 156: DICOM connection mode  Click Save if changes have been made for each entry.		☐
■ If all entries have been checked, click Finish.		☐
■ Under DICOM, select Print Devices and press under Host Properties >.		☐

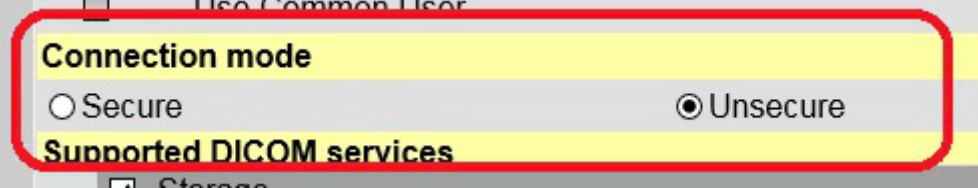ACTION	MRAWP
On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unsecure. <i>Fig. 157: DICOM connection mode</i>  Click Save if changes have been made for each entry.	☐
If all entries have been checked, click Finish.	☐
Under DICOM, select HIS/RIS Nodes and press under Host Properties >.	☐
On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unsecure. <i>Fig. 158: DICOM connection mode</i>  Click Save if changes have been made for each entry.	☐
If all entries have been checked, click Finish.	☐
Click Home and confirm the requested restart of application with OK.	☐
Wait until the restart is performed and the syngo main UI is up and running.	☐



If a new Network Node, Print device or HIS/RIS Node is configured, be aware that the default setting for the Connection Mode is secure. Please change it to Unsecure.

3.2.1.12 MRAWP (MRC) Function Test of the DICOM Nodes

ACTION	MRAWP
Send an image from the Patient Browser to all connected DICOM nodes. The images must be sent successfully	☐

ACTION	MRAWP
Send a film sheet to all connected DICOM cameras. The film sheet must be sent successfully to the DICOM camera.	☐
Retrieve a test patient from the RIS/HIS system, if available. The patient should be available on the system.	☐



A DICOM connection can now be made secure (encrypted). This is a new configuration parameter. See (→ New feature for DICOM node configuration / Page 233)



The new DICOM Node Configuration Option "Supports Presentation State" must not be enabled for any DICOM storage node. The NUMARIS/4 VE11 does not currently support the "Presentation State" for created MR images.


3.2.1.13 MRAWP (MRC) System Management




Due to the CSSI process, system management cannot be activated at the factory. System management has to be activated at the customer site.

Installation and Activation Steps

ACTION	MRAWP
Select Options --> Service -> Local Service.	☐
Enter the valid service key.	☐
In the local service UI, select Configuration	☐
Check that the settings for service region (Europe, USA, or Asia) under Local Host Site Info are correct. Correct the settings, if necessary, and save the changes. This entry assigns the IP addresses for the RDIAG server and the System Management server.	☐
Under the "Local Host" Routing, type in the IP address of the router or gateway used for SRS. In some cases, this entry is the default gateway that is used for the TCP/IP LAN configuration. See log file <code>C:\MedCom\MriSiteData\ExtBaRe\routing.txt</code> for the routing entry. For the installation of System Management, the correct entry under routing is essential. If this routing entry is not set (empty), you will get an error popup when entering the System Management configuration.	☐

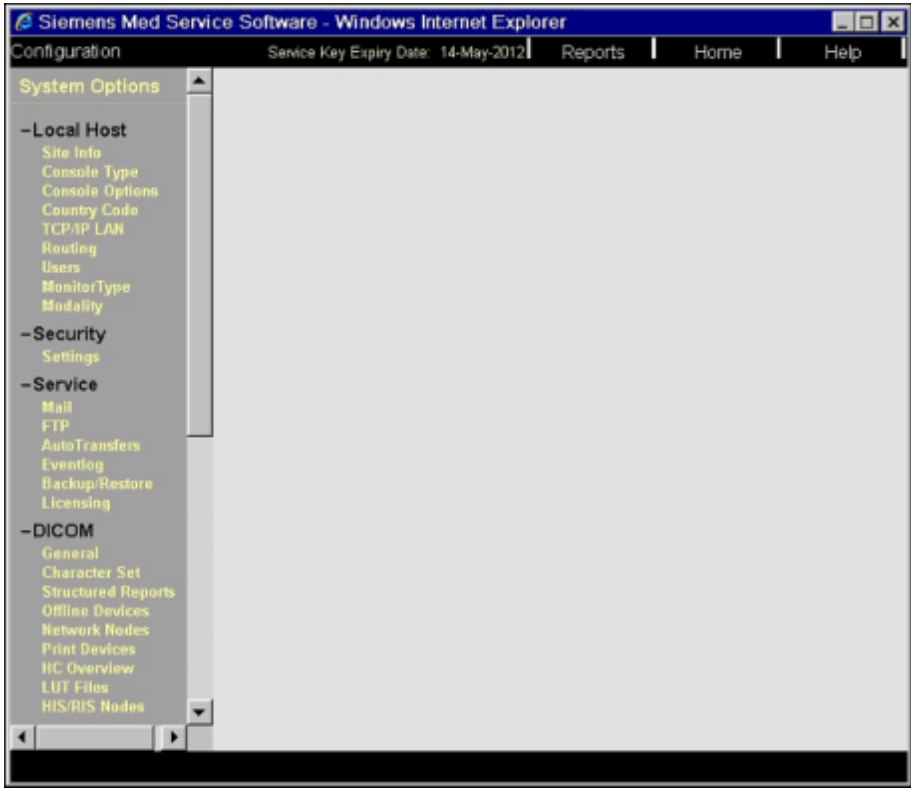ACTION	MRAWP
■ Select System Options. Fig. 159: Configuration 	☐

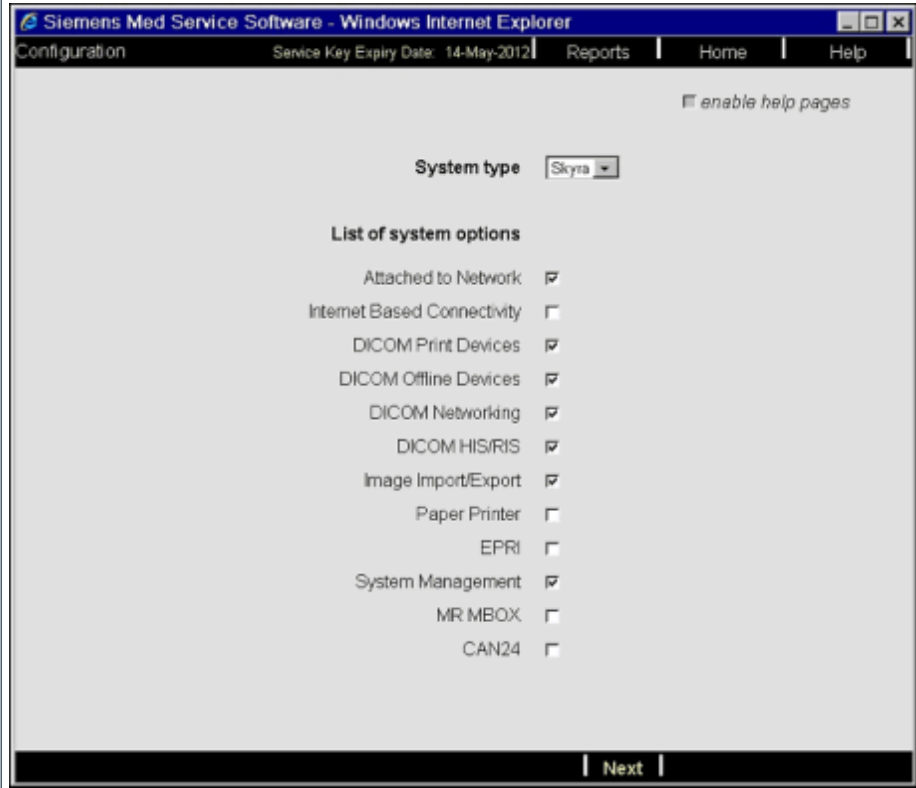
ACTION	MRAWP
■ Under "List of System Options", click the System Management checkbox Fig. 160: System options The screenshot shows a web browser window titled "Siemens Med Service Software - Windows Internet Explorer". The address bar shows "Configuration". The page has a header with "Service Key Expiry Date: 14-May-2012" and navigation links for "Reports", "Home", and "Help". There is a checkbox labeled "enable help pages". Below this is a "System type" dropdown menu set to "Skyra". The main section is titled "List of system options" and contains a list of features with checkboxes: "Attached to Network" (checked), "Internet Based Connectivity" (unchecked), "DICOM Print Devices" (checked), "DICOM Offline Devices" (checked), "DICOM Networking" (checked), "DICOM HIS/RIS" (checked), "Image Import/Export" (checked), "Paper Printer" (unchecked), "EPRI" (unchecked), "System Management" (checked), "MR MBOX" (unchecked), and "CAN24" (unchecked). At the bottom of the list is a "Next" button.	☐
■ Select Next.	☐

ACTION

MRAWP

- On the System Management Master configuration page


Fig. 161: System Management master page

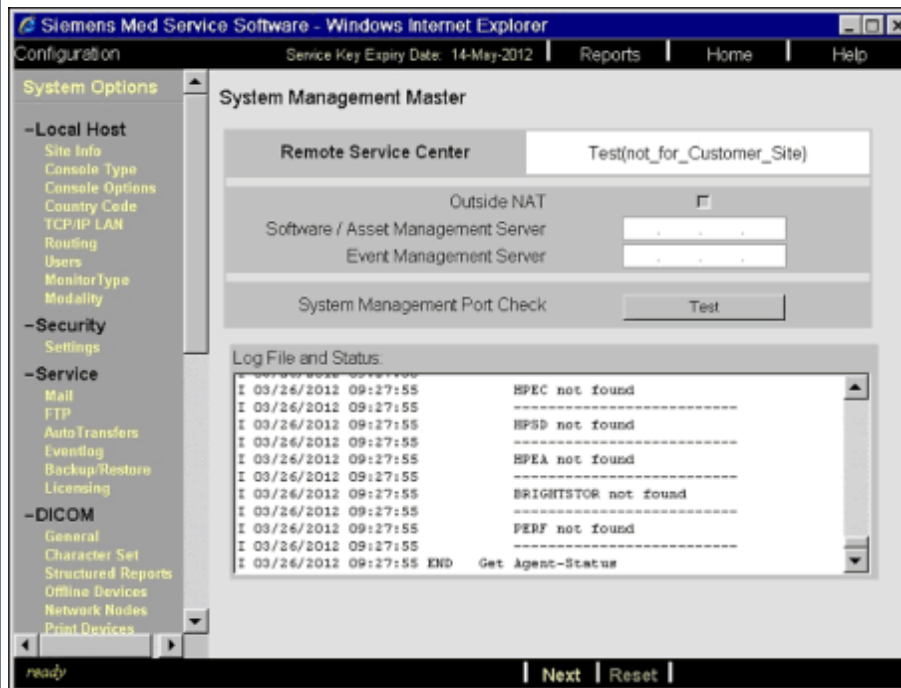

See the report file

C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html
for the Outside NAT configuration.

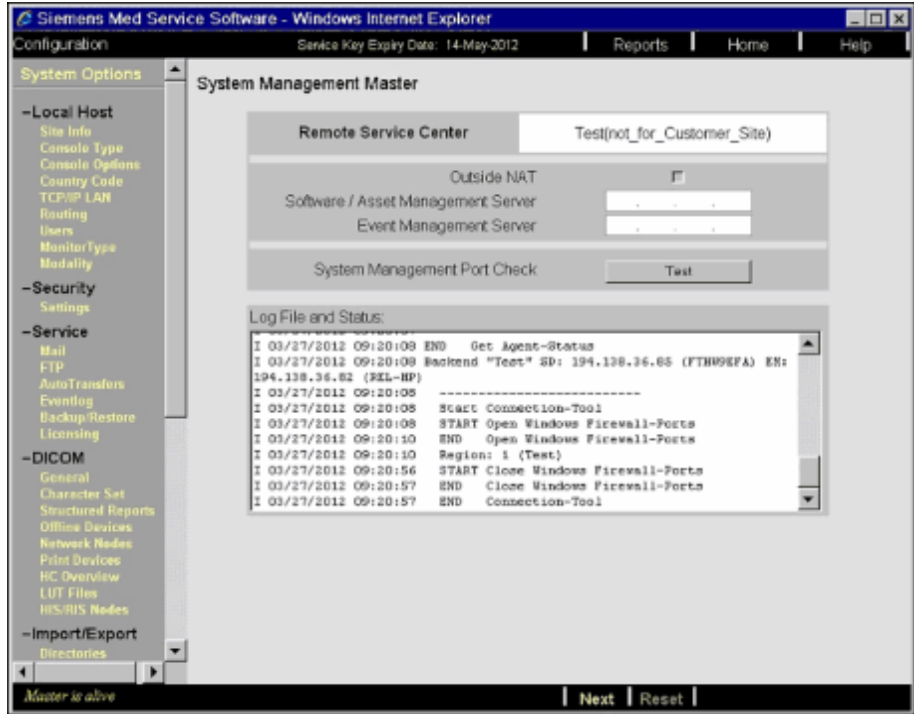

ACTION	MRAWP
■ Use the Test button to check whether the System Management server can be accessed. Fig. 162: System Management master can be reached  The log file for the test can be found under the directory: C:\MedCom\service\mwtools\CCT\Logs	☐
In the event of problems with the connection test, contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	
■ Select Next.	☐


ACTION

MRAWP

- Under **System Mgmt. Agent Controls**, set the following configuration:
Mark the checkboxes "Go Active" and "Go Online".
(The other checkboxes have already been marked).

☐

Fig. 163: System Management activation

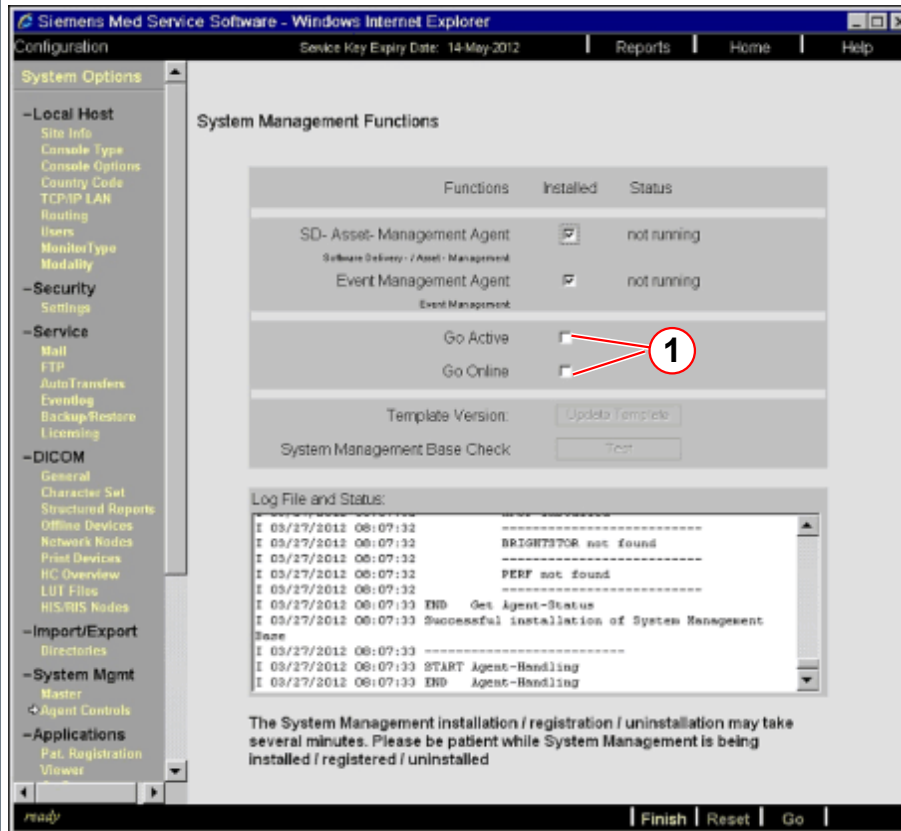


(1) Mark both checkboxes

In the Log File and Status area, the current status of System Management activation/installation is displayed.

- Click **Go**.

☐

Note: The installation/registration will take approximately 10 to 15 minutes.

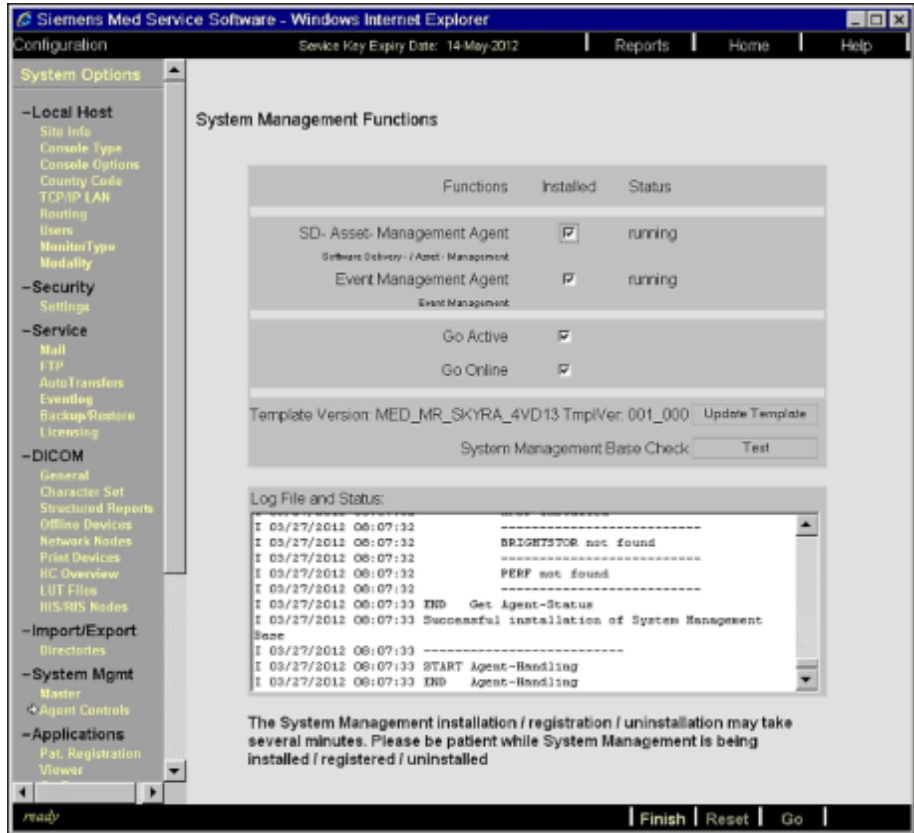ACTION	MRAWP
■ If activation is successful, the following is displayed: Fig. 164: System Management activation OK  Select Finish to continue	☐
■ Click Finish and close the local service UI. In the event of problems (installation/activation fails a second time), contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	☐



If there are problems during installation of the Event Management Agent, see (→ System management installation workaround / Page 234)

3.2.1.14 MRAWP (MRC) Restoring Changed Customer Measurement Programs from VD13/VE11A/VE11B



The following step restores all changed customer protocols. Only these protocols will be converted after the reboot.

ACTION	MRAWP
Insert the Master Backup CD labeled Master Backup VD13 or VE11A or VE11B MRAWP (actual date) created - see (→ Performing a Master Backup on the MRAWP (MRC), for changed Customer Protocols / Page 118) - into the CD/DVD-ROM drive R:	☐
Select Options --> Service --> Local Service.	☐
Enter the valid service key and click OK.	☐
Select Backup&Restore in the service platform.	☐
Select under Command Restore	☐
Select under Drives [-R-] If the drive letter is not correct, please perform (→ Drive Letter Assignment for the CD/DVD ROM Drive and/or for the CD/DVD-R Drive / MR00-000.814.55).	☐
Select under Archive the <code>MasterBackup</code> and under Groups the entry <code>N4Custom-er.</code>	☐
Press Go to start the restore. Wait until the restore has finished.	☐
Remove the Master Backup CD labeled Master Backup VD13 or VE11A or VE11B MRAWP (actual date) from the CD/DVD-ROM drive R:	☐
Click Home and confirm the requested reboot message popup with OK.	☐

3.2.1.16 MRAWP (MRC) Manual Restore of Voice Commands



Due to a bug in the backup program, it is recommended to restore the voice command files from an external medium in the new software version.



Please check the external medium first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.

ACTION	MRAWP
■ Enable the advanced user - Run the Advanced User program from the Start menu (press <CTRL> + <ESC> to access the Start menu) - From the Start menu, select All Programs -> Advanced User. - In the popup displayed, enter the advanced user password and press the OK button - Confirm Advanced User enabled popup with OK.	☐
■ Connect the USB medium containing the backup of the voice command files.	☐
■ Start Windows Explorer: - Press <CTRL> + <ESC> to access the Start menu - In the Start menu select All Programs ⇒ Accessories ⇒ Windows Explorer Windows Explorer starts.	☐
■ Navigate to the USB medium	☐
■ Copy the directory PIA and its contents from the USB medium to the directory C:\MedCom\MriSiteData	☐
■ Confirm folder replacement popup with Yes	☐
■ In the Copy File popup mark the checkbox "Do this for the next ... conflicts" and select Copy and Replace. The voice command files are now copied into the directory C:\MedCom\MriSiteData\PIA	☐
■ Close Windows Explorer.	☐
■ Remove the USB medium.	☐

3.2.1.17 MRAWP (MRC) Performing a Win 7 Image Backup on the MRAWP (MRC)



Due to the new user and password management, the previous passwords will no longer be usable. Also the user administrator will not exist. To set the passwords for the users (meduser and syadmin), see (→ New Password Handling in syngo / Page 227).

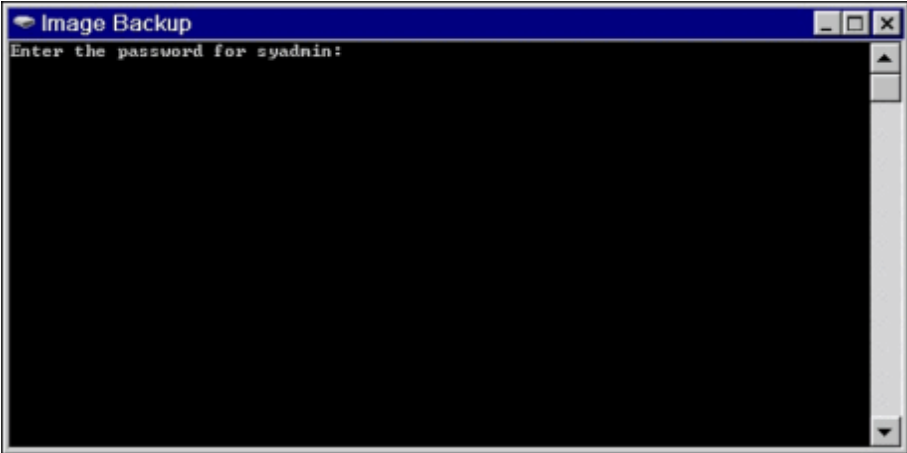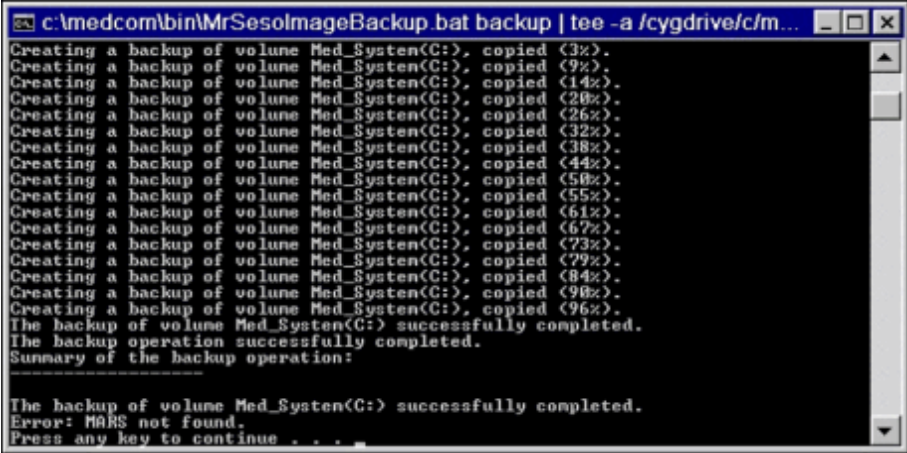ACTION	MRAWP
Create the Win 7 image backup of the system disk - Enable advanced user. - Start -> All Programs - Select Image Backup The syadmin password is requested first: Fig. 167: Password request popup  Type in the administrator password (the input is not displayed) and press ↵ Now the Windows 7 image is automatically created from the system disk, and a copy is transferred to MARS. Wait until the procedure is finished. This takes about 10 minutes. Fig. 168: Win7 image Backup  Close the Image Backup window by pressing <SPACE>	☐

3.2.1.18 MRAWP (MRC) Setup Auto Login Function



By default the autologin function is enabled. If this should be changed, see (→ New Password Handling in syngo / Page 227)

3.2.1.19 MRAWP (MRC) Additional Steps for TIM Docking Table (Mobile)



The following steps are necessary if the customer has more than one TIM docking table (mobile). The new firmware/loadware has to be loaded before the table can be used in the new Numaris version VE11C.

With the new firmware/loadware, the docking table cannot be used on MR Systems with a previous software version (VE11A or VE11B).

ACTION	MRAWP
Connect the other docking table to the MR system	☐
- A reboot of the MARS is required to start the update of the docking table. Perform the reboot as follows: - In the syngo main UI, select System - Control. - In the System Manager popup, select the Meas & Recon tab. - Press the Reboot button under Meas & Recon Control. - Confirm the displayed popup with Yes. MARS now boots.	☐
Wait until the docking table firmware update has finished. Observe the firmware/loadware update in the MR Scanner & Acquisition message popup. The update is finished when the Measurement system is ready for measurement.	☐




If no MRWP (MRSC) is available, the following steps can be skipped. Continue with (→ Final Check - to be performed for MRC and MRSC (if applicable) / Page 173)

3.2.1.20 MRWP (MRSC) Preliminary Steps before the Software Installation

ACTION	MRWP
- Ensure that the network cables at the back of the console are connected. MRSC: One network cable at the back of the computer.	☐
- Ensure that any device (e.g. monitor, in-room MRC, etc.) connected to the second video output of the host is disconnected. - Disconnect any additionally connected USB device (keep only the USB mouse and the USB keyboard connected).	☐
- If applicable, check whether the connected local printer (connected to the parallel port) is switched ON. - During the installation, the printer is automatically detected and the correct printer driver is installed.	☐

ACTION	MRWP
■ Switch on the MRSC.	☐
■ Log in as user "administrator".	☐




3.2.1.21 MRWP (MRSC) BIOS Update on MRWP (MRSC)

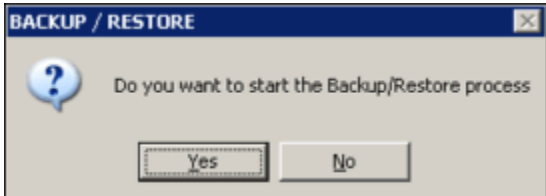
ACTION	MRWP
■ Insert the first DVD of the supplied CSS image set labeled COBR Recovery Image MR system software VE11C MR system serial no: <x> MR WP into the CD/DVD burner drive (S:\).	☐
■ Insert the supplied CD labeled with "MR-VD BIOS CD VE11.. " into the host CD/DVD-ROM Drive" (R:\)	☐
■ Invoke a command shell from the Start menu (Start ⇒ All Programs⇒ Accessories ⇒ Command Prompt)	☐
■ In the command shell, type in: S:\Start.bat	☐
■ Confirm the Installation popup with "Yes" <i>Fig. 169: Backup/Restore popup</i> 	☐
■ The running script performs a reboot at the end. This may take some time. Please wait until this is finished.	☐
■ The BIOS CD will be booted automatically.	☐
■ The BIOS update runs automatically (no user action is necessary) and ends when the following message is displayed: <pre> BIOS Installation Done - System Ready for Software Installation Please remove BIOS CD ! Please Turn Off Power or Reset System ! C:\> </pre>	☐
■ Remove the supplied BIOS CD from the CD/DVD drive (R:\)	☐




ACTION	MRWP
Remove the first DVD of the supplied CSS image set labeled COBR Recovery Image MR system software VE11C MR system serial no: <x> MR WP from the CD/DVD burner drive (S:\).	☐
Insert the first DVD of the supplied CSS image set labeled COBR Recovery Image MR system software VE11C MR system serial no: <x> MR WP into the CD/DVD-ROM drive (R:\).	☐
Press key combination <CTRL> + <ALT> + (on a German keyboard use <Strg> + <Alt> + <Entf>) to reboot and continue with "CSSI installation MRC".	☐




3.2.1.22 MRWP (MRSC) CSSI Installation

ACTION	MRWP
A reboot is performed and the COBR tool boots.	☐
Note: If you get the following error message "Operating system not found", please remove the first DVD of the supplied CSS image set from the CD/DVD drive (R:\), insert it into the CD/DVD burner drive (S:\) and perform a reboot.	

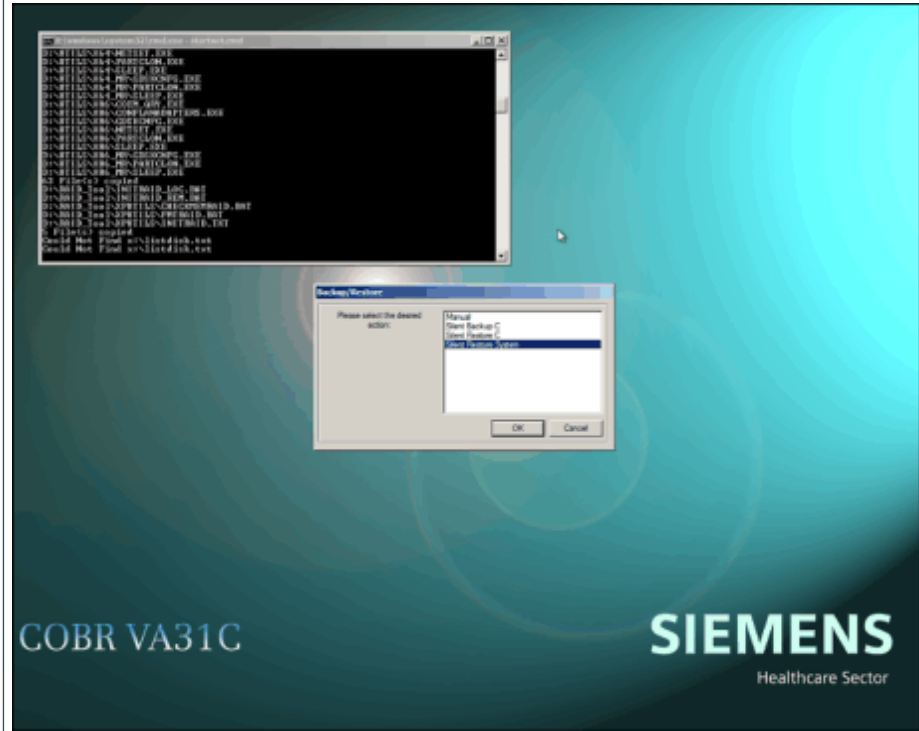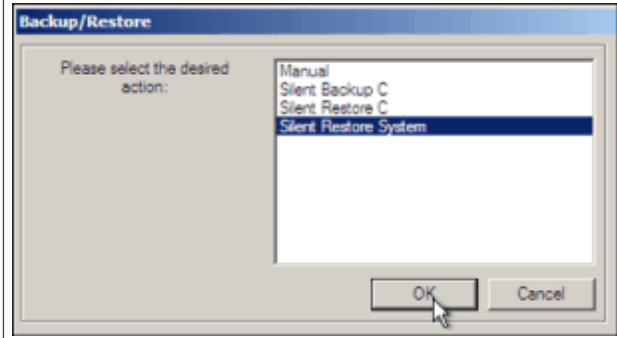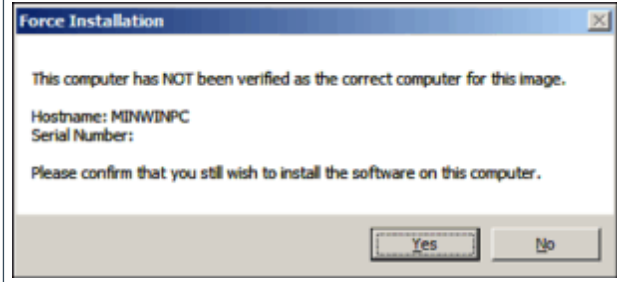ACTION	MRWP
Wait until the COBR tool window is displayed. Fig. 170: COBR tool window 	☐
In the Backup/Restore popup, select <code>Silent Restore System</code> and press OK Fig. 171: Backup/Restore popup Silent Restore System 	☐
The following popup is displayed (example is shown below) Fig. 172: Force Installation popup  Confirm this popup with Yes.	☐

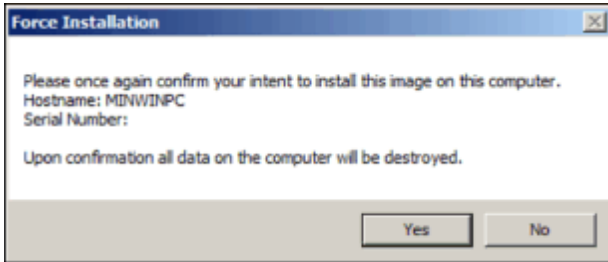ACTION	MRWP
Confirm the Force Installation popup (example is shown below) with Yes. Fig. 173: Force Installation popup  This Force Installation popup is displayed twice. Confirm it each time with Yes.	<input data-bbox="1394 271 1433 309" type="checkbox"/>
The system restore process starts now. The following shows an example of the screen display. Fig. 174: COBR restore process start 	<input data-bbox="1394 680 1433 719" type="checkbox"/>

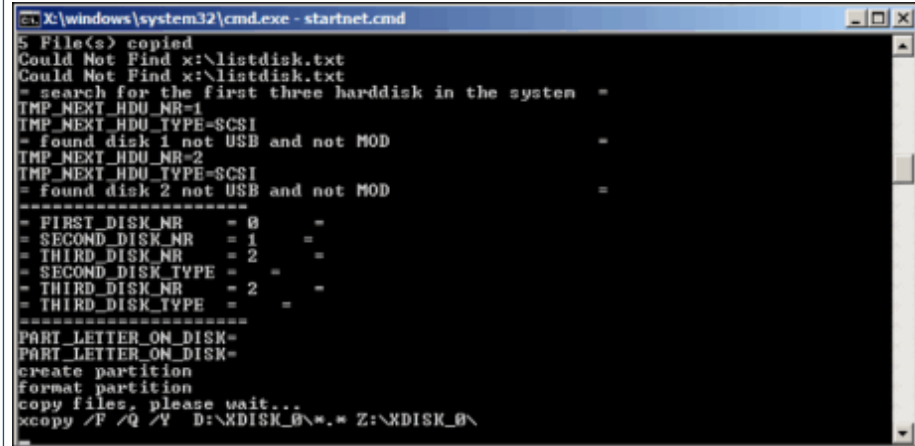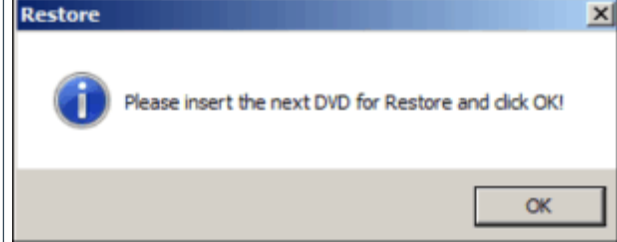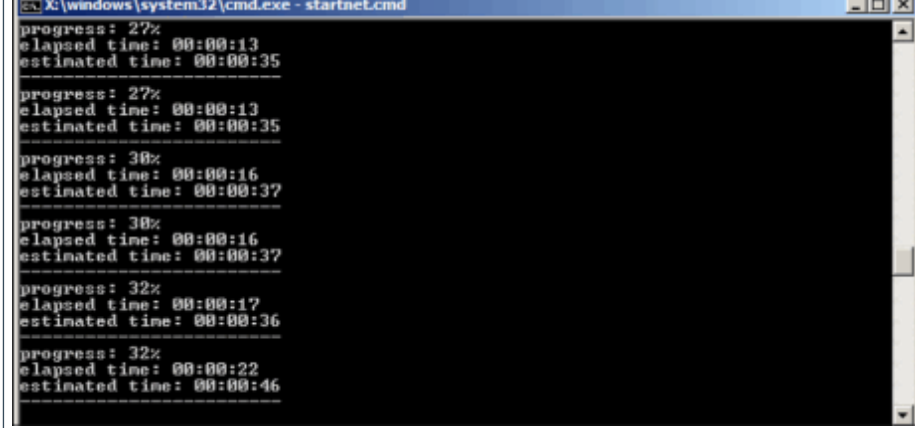ACTION	MRWP
First all image files from the DVD for the system disk (disk1) are copied from the DVD to the hard disk. Fig. 175: Copy image files  After the files are copied, the DVD is automatically ejected.	☐
Remove the ejected medium from the CD/DVD-ROM drive, insert the next DVD from the system image DVD set, and close the drive.	☐
Confirm the request popup with OK. Fig. 176: Next DVD request popup 	☐
The copy procedure continues. The next DVD will be requested automatically until all required data has been copied. The copy and the confirmation of the popup (→ Fig. 176 Page 158) have to be repeated for all available CSSI DVDs.	
Now the image restore starts. Fig. 177: Image restore progress output 	☐

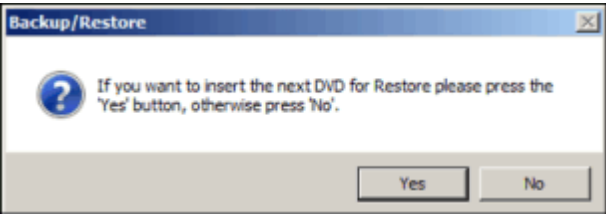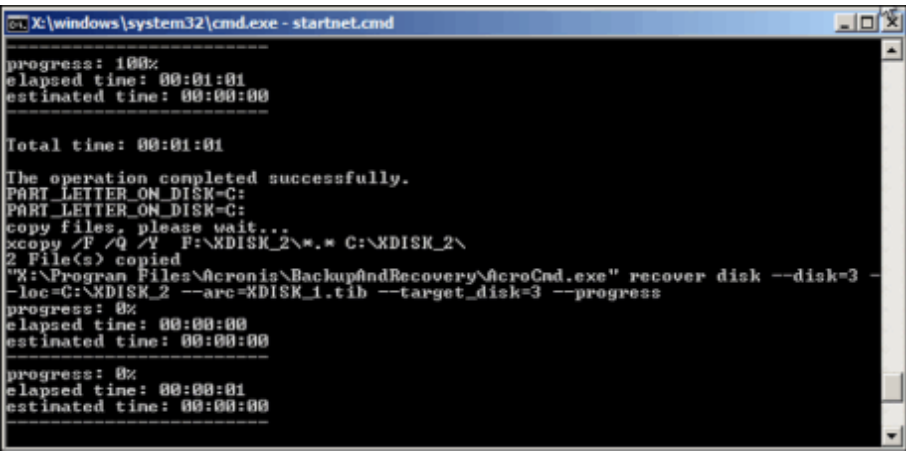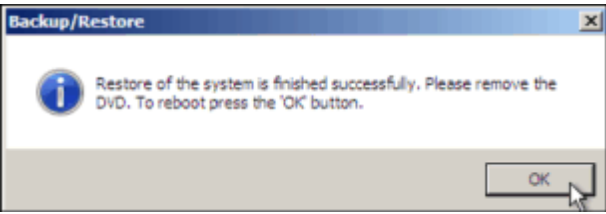
ACTION	MRWP
Wait until the following popup appears. <i>Fig. 178: Next DVD request popup</i> 	☐
Confirm the "Next DVD request popup" shown above with "No".	☐
Now the last disk is restored from the image. <i>Fig. 179: Last disk restore</i> 	☐
When all data has been successfully restored, the following popup is displayed. <i>Fig. 180: Backup/Restore successfully finished popup</i> 	☐
Remove the ejected medium, close the drive, and confirm this popup with OK.	☐
The system is rebooted.	☐
Note: The following behavior is new in the VE11C CSSI update.	
After the reboot, the autologin fails with an error message. Please wait until a further reboot is done (takes about 1 minute) The next autologin will be successful. Note: Please don't click or hit any key - please just wait.	☐
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 5 minutes. Note: Don't close or confirm any displayed popups.	☐

ACTION

MRWP

During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:


Fig. 181: Displayed popup

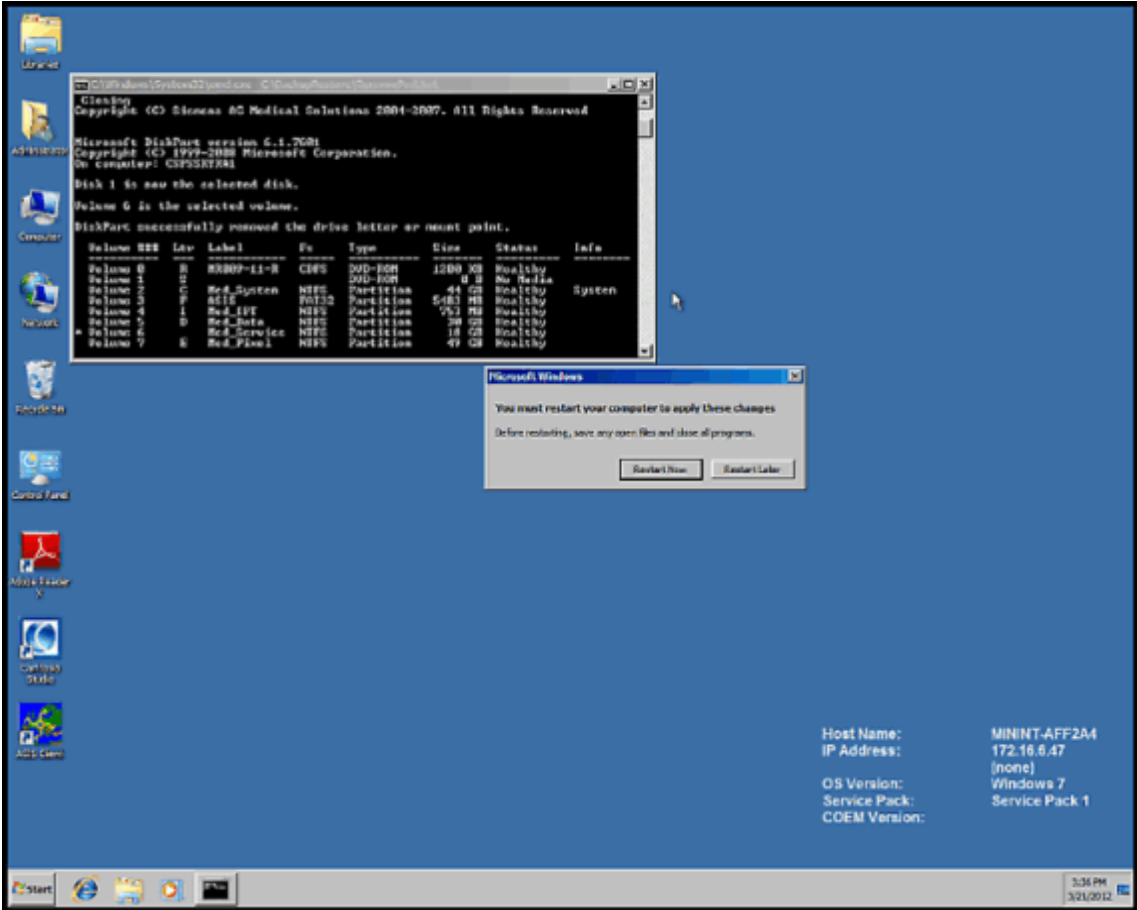



The CSSI Post popup is also displayed (see Fig. 60 "CSSI Post Actions popup")

Note: Don't accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don't touch or accept it. Keep it displayed on the screen. It will be closed automatically.

Fig. 182: Windows popup

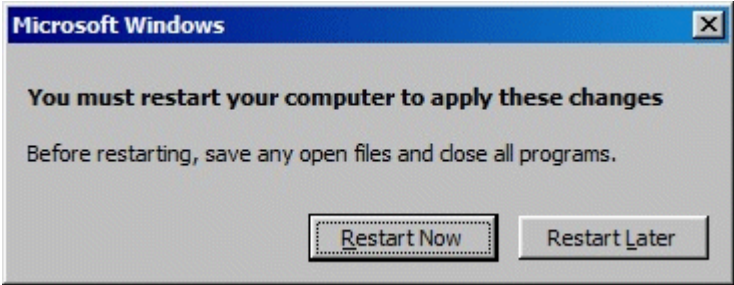

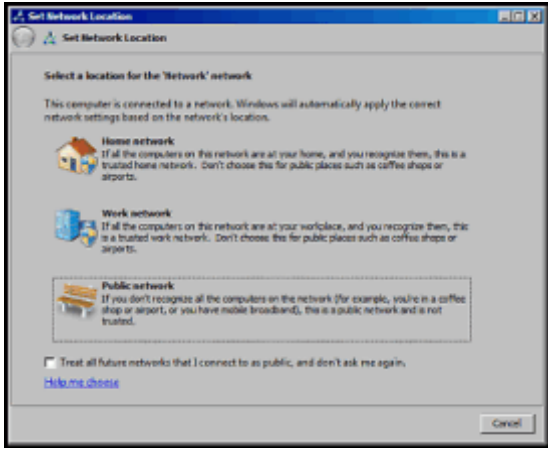ACTION	MRWP
Note: The following popup can be ignored too, if it is displayed. Fig. 183: Network location popup 	
■ The system is rebooted.	☐
■ An autologin as user "syadmin" is performed and the following CSSI post popup is displayed.	☐

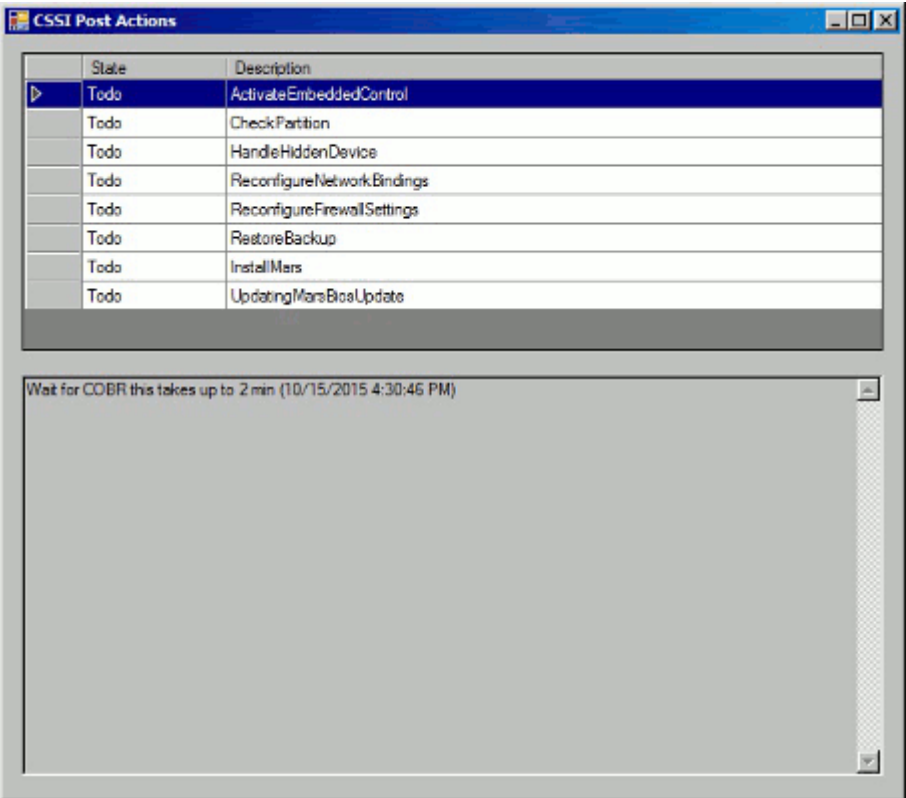ACTION	MRWP
The CSSI post popup shows the status and progress of the different post CSSI update steps. <i>Fig. 184: CSSI Post Actions popup</i>  Please wait until the steps are finished. During these post CSSI update steps the host will be rebooted automatically. Note: No user interaction is necessary.	☐
The system performs an autologin as meduser after all post CSSI update steps are done.	☐
Wait until the syngo main UI is up and running.	☐

3.2.1.23 MRWP (MRSC) Checking the Patient Browser Connection with Shared Database



Ensure that the satcon with remote database configuration is not configured as a network node to the Main console. Otherwise problems will occur.

ACTION	MRWP
Use the Patient Browser to check whether the MRC database can be reached.	☐

3.2.1.24 MRWP (MRSC) Network Printer Installation on the MRSC

ACTION	MRWP
Install the network paper printer according to the software installation manual (this step is not part of the dynamic data) if previously installed on the MRAWP or MRWP.	☐
If applicable, set the installed printer as the default printer. Select Start -> Settings -> Printers and Faxes Select the printer and click the right mouse button. In the menu displayed, select the "Default printer" entry.	☐

3.2.1.25 MRWP (MRSC) Configuration checks

ACTION	MRWP
Select Options ⇒ Service ⇒ Local Service.	☐
Enter the valid service key and click OK.	☐
Select Configuration in the service platform.	☐

ACTION	MRWP
■ Select Monitor Type under Local Host and check the settings as follows: - Monitor Type: Color and Siemens (True Color) - Multi-Monitor settings Only "Start syngo in single monitor" is marked Fig. 185: Monitor settings SESO 	☐
■ Click Save if changes have been made and Finish	☐

ACTION	MRWP
On the next page, under General Node Properties in the drop-down menu Logical Name, select each available entry. For each entry check the Connection Mode. This must be set to Unsecure. Fig. 188: DICOM connection mode  Click Save if changes have been made for each entry.	☐
If all entries have been checked, click Finish.	☐
Click Home and confirm the requested restart of application with OK.	☐
Wait until the restart is performed and the syngo main UI is up and running.	☐



If a new Network Node, Print device or HIS/RIS Node is configured, be aware that the default setting for the Connection Mode is secure. Please change it to Unsecure.

3.2.1.26 MRWP (MRSC) Function Test of the DICOM Nodes

ACTION	MRWP
Send an image from the Patient Browser to all connected DICOM nodes. The images must be sent successfully	☐
Send a film sheet to all connected DICOM cameras. The film sheet must be sent successfully to the DICOM camera.	☐



A DICOM connection can now be made secure (encrypted). This is a new configuration parameter. See (→ New feature for DICOM node configuration / Page 233)




The new DICOM Node Configuration Option "Supports Presentation State" must not be enabled for any DICOM storage node. The NUMARIS/4 VE11 does not currently support the "Presentation State" for created MR images.

3.2.1.27 MRWP (MRSC) System Management



Due to the CSSI process, system management cannot be activated at the factory. System management has to be activated at the customer site.

Installation and Activation Steps

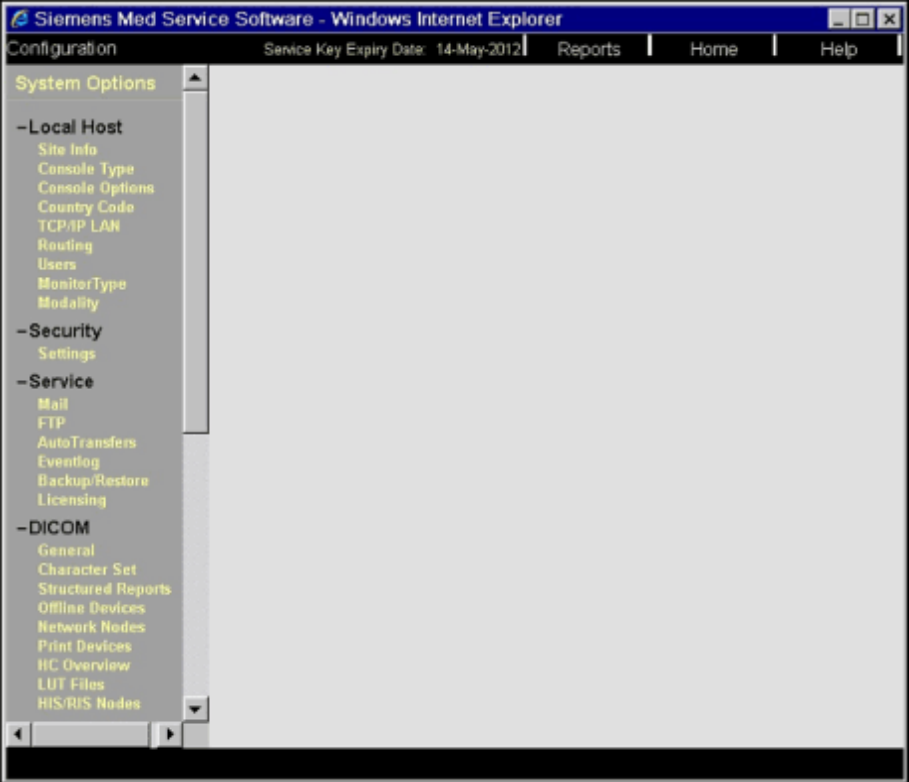ACTION	MRWP
Select Options --> Service -> Local Service.	☐
Enter the valid service key.	☐
In the local service UI, select Configuration	☐
Check that the settings for service region (Europe, USA, or Asia) under Local Host Site Info are correct. Correct the settings, if necessary, and save the changes. This entry assigns the IP addresses for the RDIAG server and the System Management server.	☐
Under the "Local Host" Routing, type in the IP address of the router or gateway used for SRS. In some cases, this entry is the default gateway that is used for the TCP/IP LAN configuration. See log file C:\MedCom\MriSiteData\ExtBaRe\routing.txt for the routing entry. For the installation of System Management, the correct entry under routing is essential. If this routing entry is not set (empty), you will get an error popup when entering the System Management configuration.	☐
Select System Options. Fig. 189: Configuration 	☐

ACTION	MRWP
■ Under "List of System Options", click the System Management checkbox Fig. 190: System options The screenshot shows a web browser window titled "Siemens Med Service Software - Windows Internet Explorer". The address bar shows "Configuration". The page has a navigation bar with "Service Key Expiry Date: 14-May-2012", "Reports", "Home", and "Help". Below the navigation bar, there is a checkbox labeled "enable help pages". The main content area is titled "System type" with a dropdown menu showing "Skyra". Below this is a section titled "List of system options" containing a list of options with checkboxes: "Attached to Network" (checked), "Internet Based Connectivity" (unchecked), "DICOM Print Devices" (checked), "DICOM Offline Devices" (checked), "DICOM Networking" (checked), "DICOM HIS/RIS" (checked), "Image Import/Export" (checked), "Paper Printer" (unchecked), "EPRI" (unchecked), "System Management" (checked), "MR MBOX" (unchecked), and "CAN24" (unchecked). At the bottom of the window, there is a "Next" button.	☐
■ Select Next.	☐


ACTION

MRWP

- On the System Management Master configuration page




Fig. 191: System Management master page

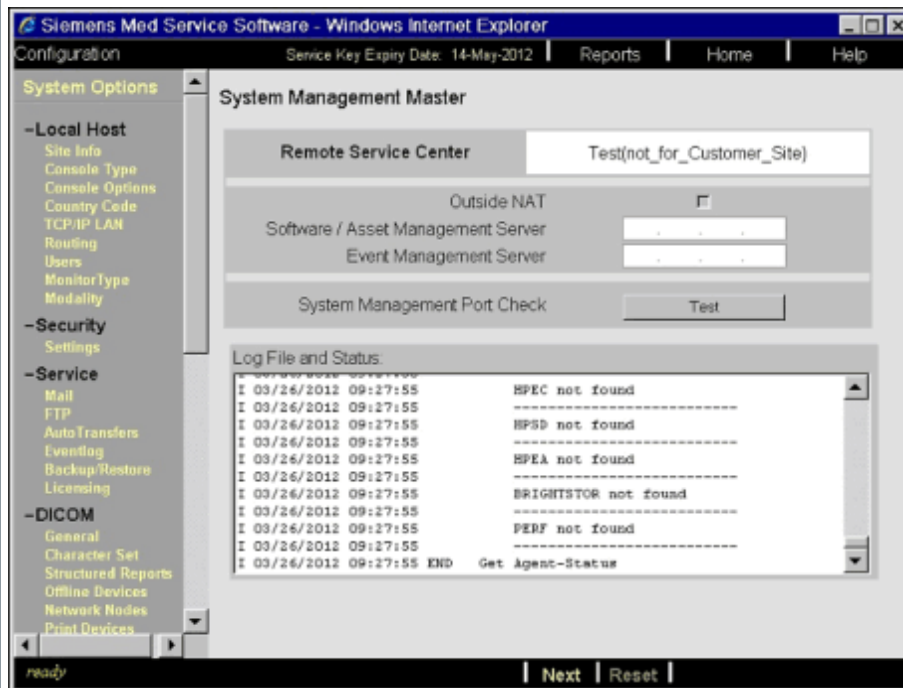

See the report file

C:\MedCom\MriSiteData\BaRe_HW_Info\MRC\report.html
for the Outside NAT configuration.

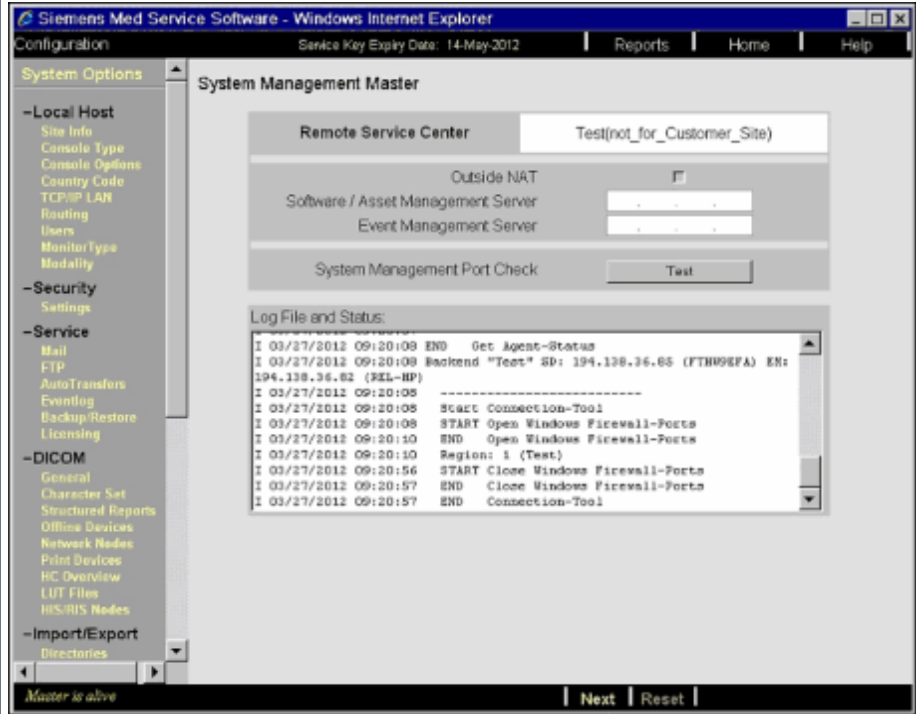ACTION	MRWP
■ Use the Test button to check whether the System Management server can be accessed. Fig. 192: System Management master can be reached  The log file for the test can be found under the directory: C:\MedCom\service\mwtools\CCT\Logs	☐
In the event of problems with the connection test, contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	
■ Select Next.	☐


ACTION

MRWP

- Under **System Mgmt. Agent Controls**, set the following configuration:
Mark the checkboxes "Go Active" and "Go Online".
(The other checkboxes have already been marked).

☐

Fig. 193: System Management activation

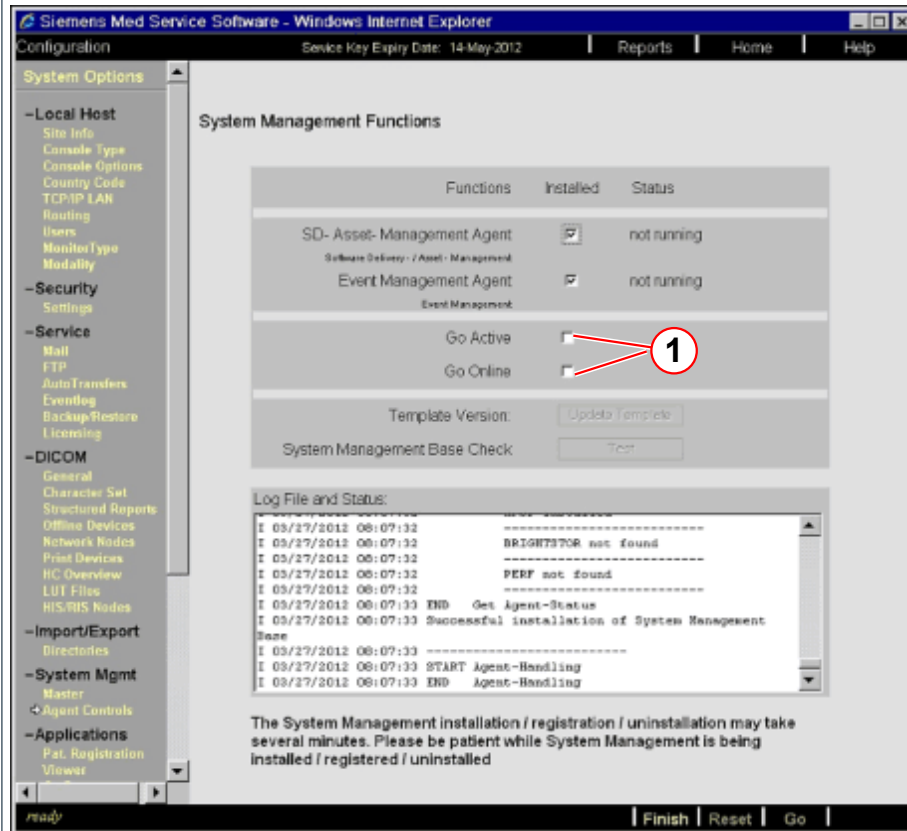


(1) Mark both checkboxes

In the Log File and Status area, the current status of System Management activation/installation is displayed.

- Click **Go**.

☐

Note: The installation/registration will take approximately 10 to 15 minutes.

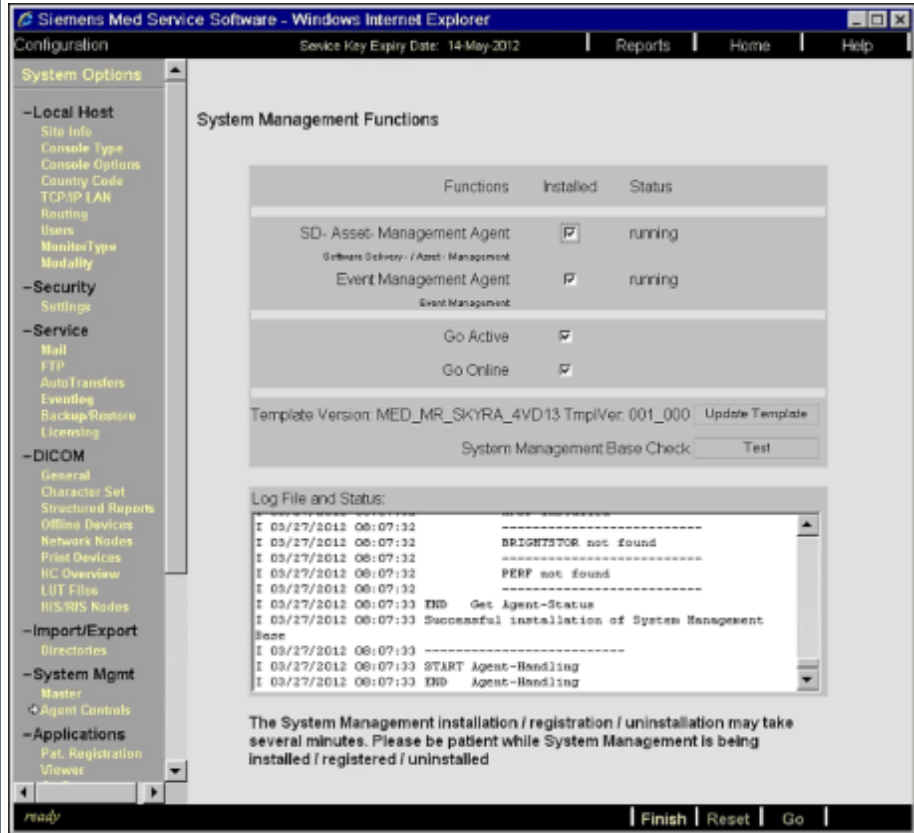ACTION	MRWP
■ If activation is successful, the following is displayed: Fig. 194: System Management activation OK  Select Finish to continue	☐
■ Click Finish and close the local service UI. In the event of problems (installation/activation fails a second time), contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.	☐



If there are problems during installation of the Event Management Agent, see (→ System management installation workaround / Page 234)

3.2.1.28 MRWP (MRSC) Setup Auto Login Function



By default the autologin function is enabled. If this should be changed, see (→ New Password Handling in syngo / Page 227)

3.2.1.29 MRWP (MRSC) Setup Advanced User Login Function



By default the Advanced User Login function is disabled. If this should be changed, see (→ Advanced User handling / Page 231)



Please contact the customer if the advanced user should be enabled, and ask the customer/admin what should be used as the advanced user password.

Also ask the customer/admin what should be used for other users. Then hand over the passwords to the customer/admin (See (→ New Password Handling in syngo / Page 227) and (→ Advanced User handling / Page 231))

3.2.1.30 MRWP (MRSC) Perform a COBR Image Backup on the MRWP (MRSC)



Please check the external medium (e.g. USB disk) first with the virus scanner of the service PC. Ensure that the external medium used is free of viruses.



A Win7 image backup with the script from the Start menu is not recommended due to the limited size of the partition D.

ACTION	MRWP
■ Create the COBR image backup of the system accordingly - Manual Backup of System Images to USB-Disk or - Manual Backup of System Images to DVD Note: Due to performance and time issues, it is recommended to use a USB disk for backup in place of DVDs. The COBR image tool must be started in the following way due to Embedded Control used in VE11C: - Login as user syadmin - Invoke a command shell from the Start menu (Start ⇒ All Programs⇒ Accessories ⇒Command Prompt) - Type in the command: <code>C:\MedCom\bin\StartDiskImageGeneration.exe ↵</code> Follow the instruction displayed in the cmd tool window.	☐

3.2.1.31 Final Check - to be performed for MRC and MRSC (if applicable)



If the customer uses a Network Time Protocol Server for synchronizing the time on the MR system, the configuration has to be done again. This configuration cannot be done in the factory due to the check to determine if the NTP server is reachable (local service UI, Site Info, step 2)

- Verify the correct installation as follows:

In the main syngo MR UI, select

Help -> Info

Check whether the following entry is visible under the Siemens software package

- N4_VE11C_LATEST_20160120

The label

- N4_VE11C_LATEST_20160120

indicates a correct N4VE11C base installation.

- Check the RDIAG configurations on the **MRAWP** (MRC) and, if applicable, on the **MRWP** (MRSC) as well (the MRSC must also be connected to RDIAG) as follows:
 - Check the host name and IP address of the configured RDIAG server (Local Service -> Configuration -> Service Mail)

	Host name	IP address
Europe, Africa	LUX09505	194.138.39.18
North/South America	SRSACC1	129.73.116.92
Asia	sgpt806x	194.138.243.178

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- In the second step of the mail configuration, check the defined mail target (for the STE, the e-mail address must be changed according to the following table)

	Logical name for FTP target	Select the mail server host	E-mail address
Europe, Africa	SRS_mail	LUX09505	MRSMC.healthcare@siemens.com
	SaveLog	LUX09505	Savelog.healthcare@siemens.com
North/South America	SRS_mail	SRSACC1	MRUPTIME.healthcare@siemens.com
	SaveLog	SRSACC1	Savelog.healthcare@siemens.com
Asia	SRS_mail	sgpt806x	MRSMC.healthcare@siemens.com
	SaveLog	sgpt806x	Savelog.healthcare@siemens.com

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- In the second step of the FTP configuration (Local Service -> Configuration -> Service FTP), check the defined FTP target

	Logical name for FTP target	FTP host selection
Europe, Africa	SRS_ftp	LUX09505

	Logical name for FTP target	FTP host selection
North/South America	SRS_ftp	SRSACC1
Asia	SRS_ftp	sgpt806x

FTP user name: `rdiagftp`

Password: `siemens`

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Test the FTP connection to the RDIAG server as follows:

Local Service -> File&Image Tools

Source: `File`

Action: `Transfer to System`

File Type: `Any file`

From source: `SRS_ftp`

If a list from the directory `/outgoing/mr` (directory on the RDIAG server) is displayed, the connection and FTP access are OK

If the test is not OK, check the SRS configuration; refer to

- Click Home
- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **EL_Sample**. The following settings must be applied:

Active checkbox: `must be marked`

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **EL_Transfer**. The following settings must be applied:

Active checkbox: `must be marked`

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **IB_Sample**. The following settings must be applied:

The "Activated" checkbox `must be marked`

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **IB_Transfer**. The following settings must be applied:

The "Activated" checkbox `must be marked`

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **edb**. The following settings must be applied:

Active checkbox: must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- Check the Autoreport setting (Local Service -> Autoreport) for the report with the name **qdb**. The following settings must be applied:

Active checkbox: must be marked

Make any necessary corrections (a detailed description is available in the Remote Connectivity chapter of the Software Installation manual).

- The **UTILDATA** auto report is used for transferring the utilization data of the system. Utilization Management is optional; for legal reasons the customer has to sign a contract in advance. Please consult your service organization before you activate System Utilization.



Prevent failure messages where the regional language setting did not match, e.g. UI language

- guidance text used within the workflow (Siemens pre-defined texts)
- language text used in error messages
- table position of the protocol after opening in a completely different range, and not useful for planning. (Due to mismatch between "." and ",")

See SKB0022441 entry and (→ Failure messages where the regional language setting did not match, e.g. UI language / MR00-000.814.55).

3.2.1.32 Customer Information

The mandatory (x days on-site) Appl. Training is being provided via MR00x/16/E.

Any additional training required has to be determined individually.

3.2.1.33 Removed Windows Feature in VE11C

Due to OS hardening, the following Windows features are removed from the VE11C.

- Windows Gadget platform
- Windows search
 - Search functionality in task menu
 - Search functionality in Windows Explorer
- XPS Viewer
- Media features
 - Windows DVD maker
 - Windows Media Center
 - Windows Media Player
- Remote Differential Compression

- Tablet PC Components
- Games
- Simple TCP Services
- TFTP Client

3.2.1.34 Final Work Steps

1. **For systems equipped with the Expert-i license:**
Follow the instructions in document MR00-000.814.28 to update existing client PCs using the Expert-i client software delivered with the CSSI upgrade.
Otherwise scrap the Expert-i client software.
2. Update the system documentation. Fill out, and if necessary, make a copy of the attached "Completion Protocol" and file it in the corresponding folder.
3. The update has to be reported by authorized personnel, either via an intranet application or by automated data transfer (only for countries connected via the update handling interface).

3.2.1.35 Changes to Previous Version

n.a.

4.1 Appendix

4.1.1 Drive Letter Assignment for the CD/DVD ROM Drive and/or for the CD/DVD-R Drive

4.1.1.1 Description:

The drive letter assignment for the CD/DVD ROM drive and/or for the CD/DVD-R drive is not correct. The following is displayed in the local service configuration

Fig. 195: Error message



Fig. 196: Error message


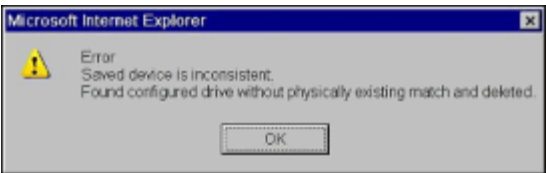


Fig. 197: Local service UI

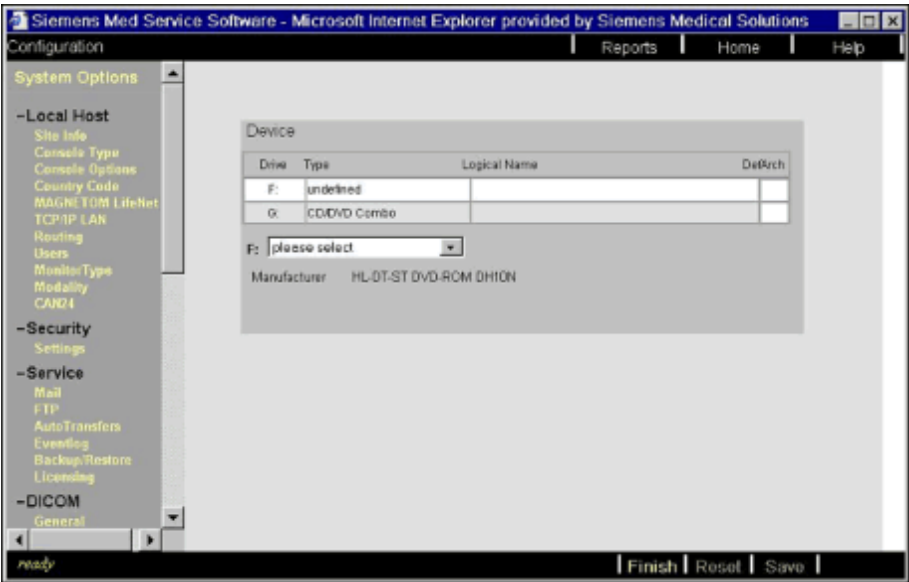


4.1.1.2 Solution:

The following procedure describes the correction for the drive letter assignment.

- Log in as the "administrator".


- Start the disk manager by selecting **START** and entering `diskmgmt.msc` ↵

Fig. 198: Start Menu

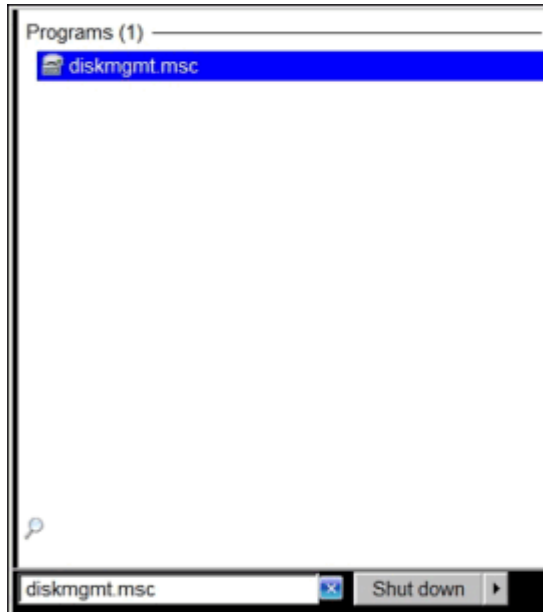

- The disk manager pops up.
- Select the drive to be (re-) assigned with a right mouse click, select '**Change drive letter and paths...**' from the pop-up menu, select '**Change**' in the next pop-up and choose the proper drive letter from the drop-down list of the next one. Click '**Yes**' in the 'Confirm' popup. (Repeat this sequence for every drive for which you have to change the drive letter).

Fig. 199: Disk management

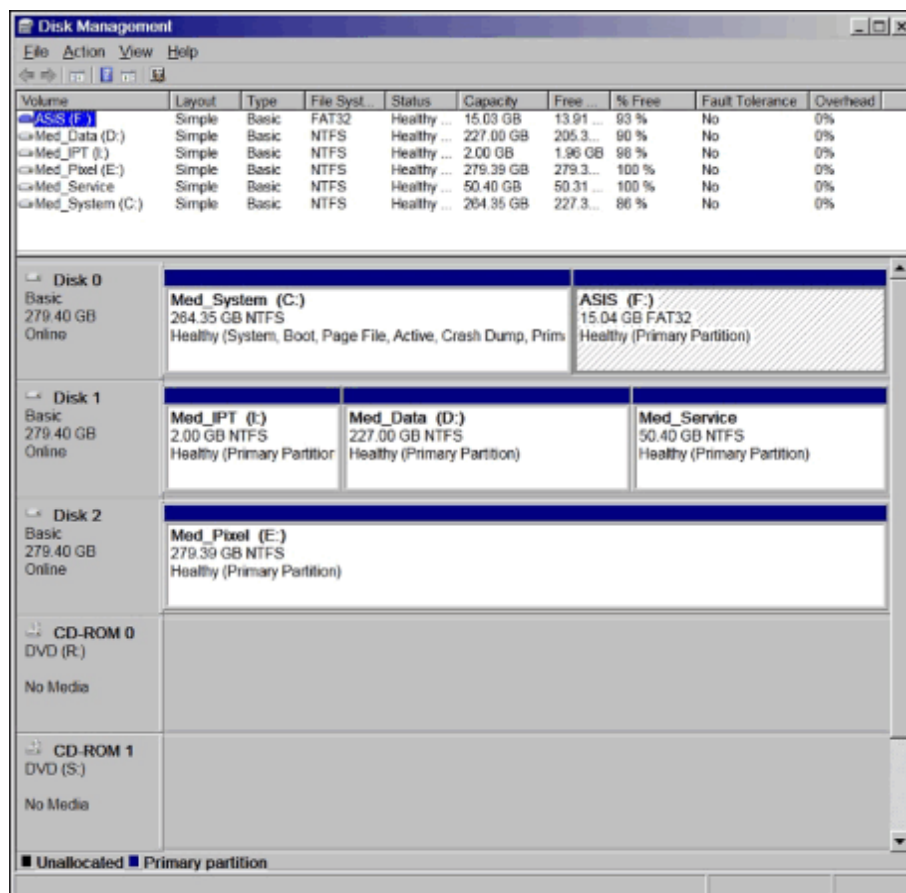

- CD/DVD-ROM drive -> **R:**
- CD/DVD-R (writer) drive-> **S:**
- Perform a reboot.
- Reconfigure Offline devices

4.1.2 Correct syngoUI Language Settings

Problem: No message displayed in the MR Scanner & Acquisition popup.


The reason is that the syngo main UI language is not correct set, wrong format is used. Please correct it see below.

Language change:

- Select in the syngo main UI: **Options** ⇒ **Configuration** ⇒ **Regional and Language Options**

Select the correct language in the Formats tab under the drop down box Format.

Fig. 200: Regional and Language Options

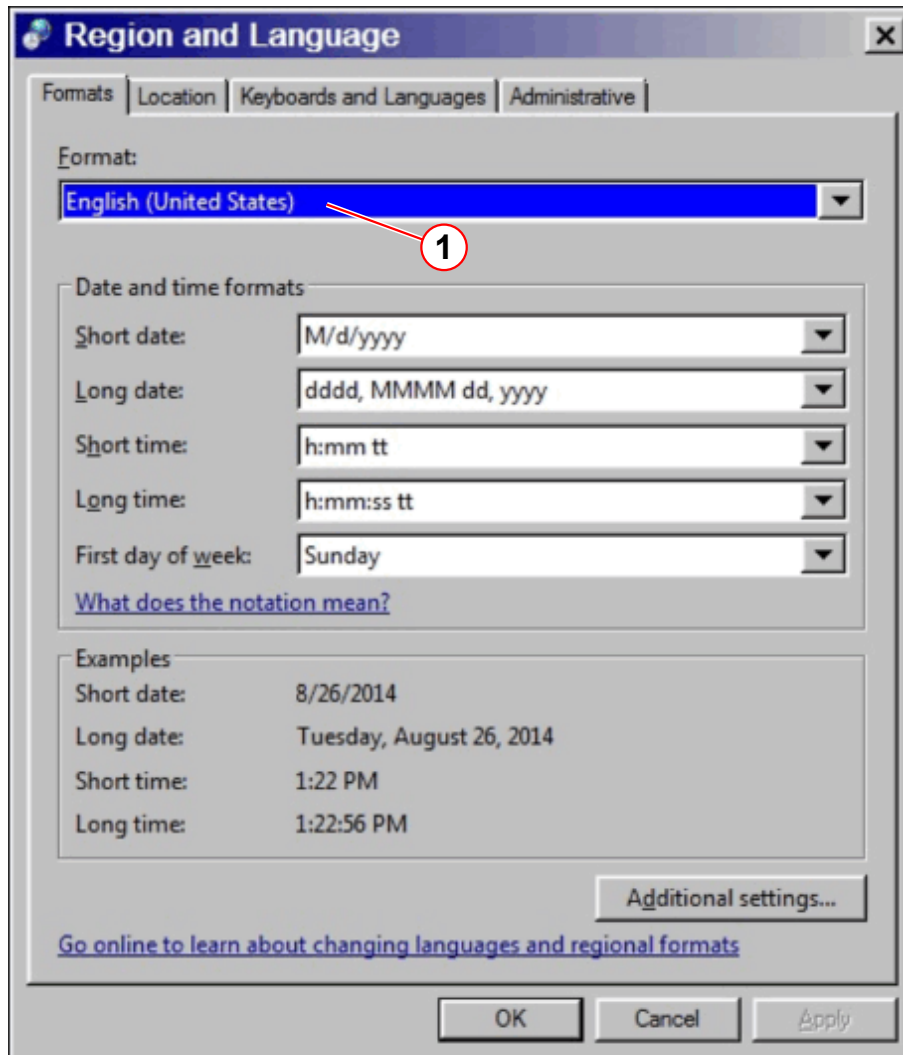

Only the following entries can be used for 1 in the graphic above:

Fig. 201:

main UI language	entry in the pull down menu
german	German (Germany)
english	English (US)
english (metric)	English (UK)
french	French (France)
spanish	Spanish (Spain)
chinese	Chinese (PRC)
japanese	Japanese

- Press OK for restart the system
- Check the scanner symbol. If the MARS isn't connected, the scanner symbol has to be crossed out. Select this symbol, and a new window with messages should appear

4.1.3 Failure messages where the regional language setting did not match, e.g. UI language

The following following failure messages can occur:

- to guidance text within the workflow (Siemens pre-defined texts)
- to language text in error messages
- table position of the protocol after opening in a total different range and not useful for planning. (Due to mismatch between "." and ",")

The following has to be checked in order to prevent these mentioned issues:

- Before doing this: Close the patient!
- Deactivate Advanced User! (start > All Programms > MR_User)

- Select from the syngo main UI Options > Configuration > Regional Settings > use the first **format** entry

Fig. 202: Example for English


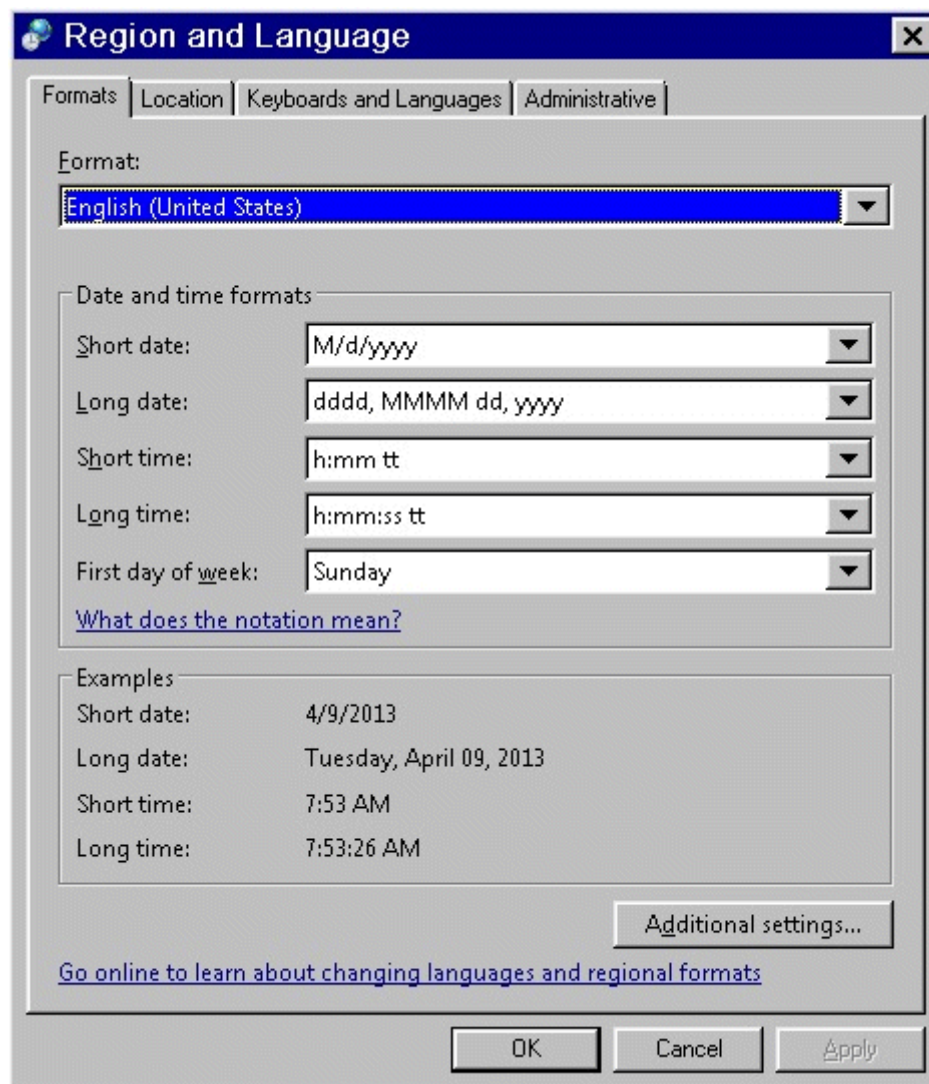


Fig. 203: Example for German



Press **Apply** and **OK**.

- After switching the language a reboot is pending. This reboot is mandatory!

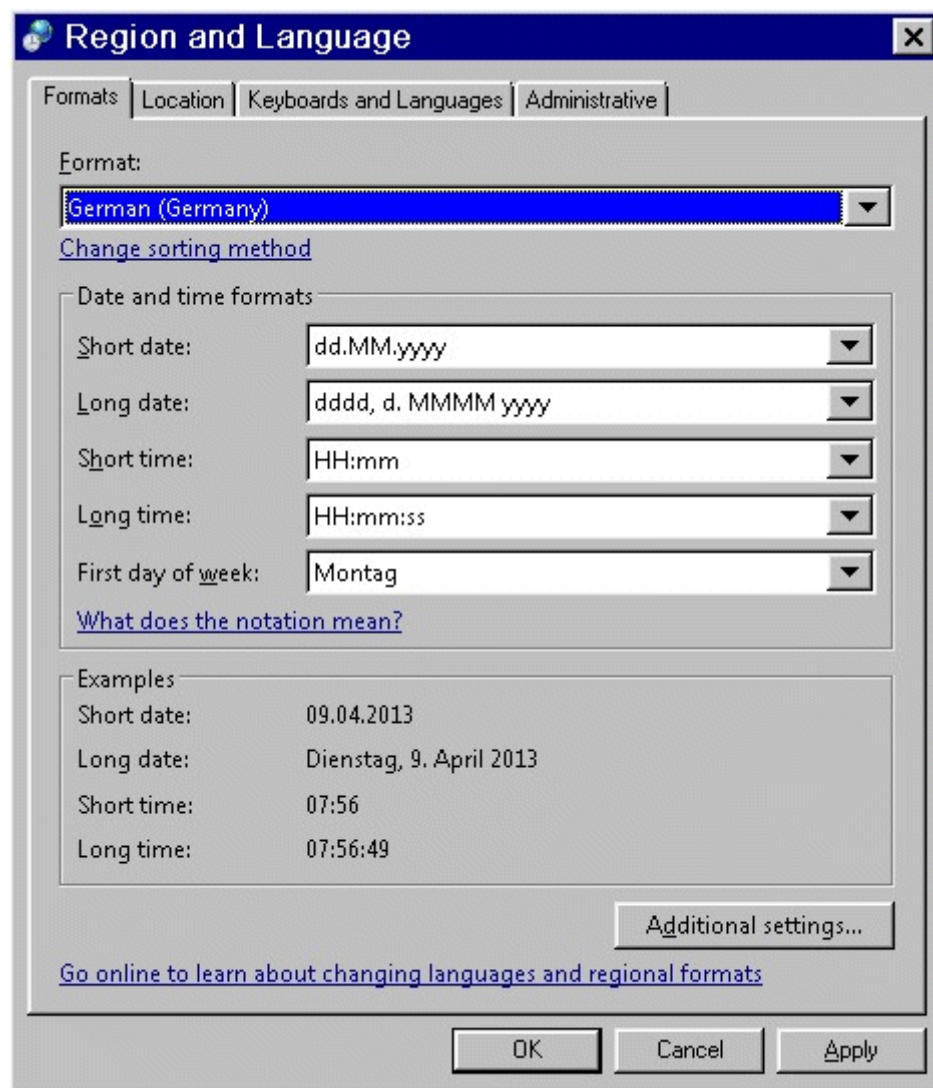

This rules are now with VE11A also valid for Japanese and Chinese language.

4.1.4 Language settings

Supported languages for meduser


- English (US)
- German (Germany)
- French (France)
- Spanish (Spain)
- Japanese

- Chinese (PRC)

4.1.4.1 Keyboard and Language Settings - Japanese


- In the Service UI, select Options -> Configuration
- Select the DICOM Character Set > in the drop-down menu DICOM default character set, select **Japanese, Katakana**

Fig. 204: Configuration - DICOM default character set

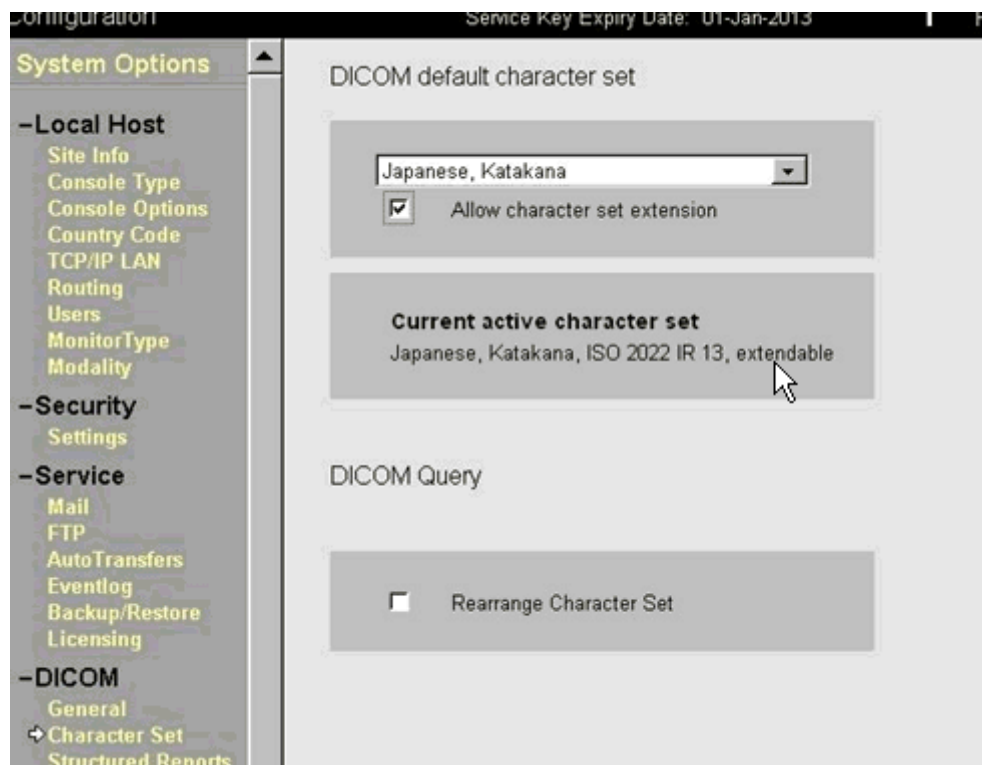

- Click **OK** to confirm the change in character set

Fig. 205: Configuration - character set data successfully saved

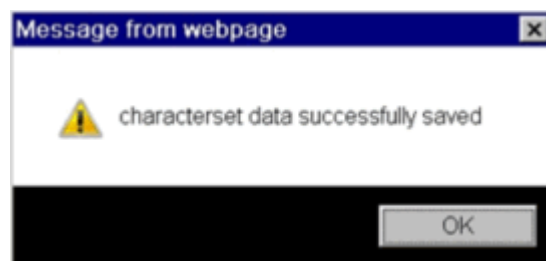


- Click **OK** to restart the application

Fig. 206: Configuration - a restart of the application SW is necessary to make changes valid




- Select **Options > Configuration**
- Select **Regional and Language Options**.

Fig. 207: Configuration panel

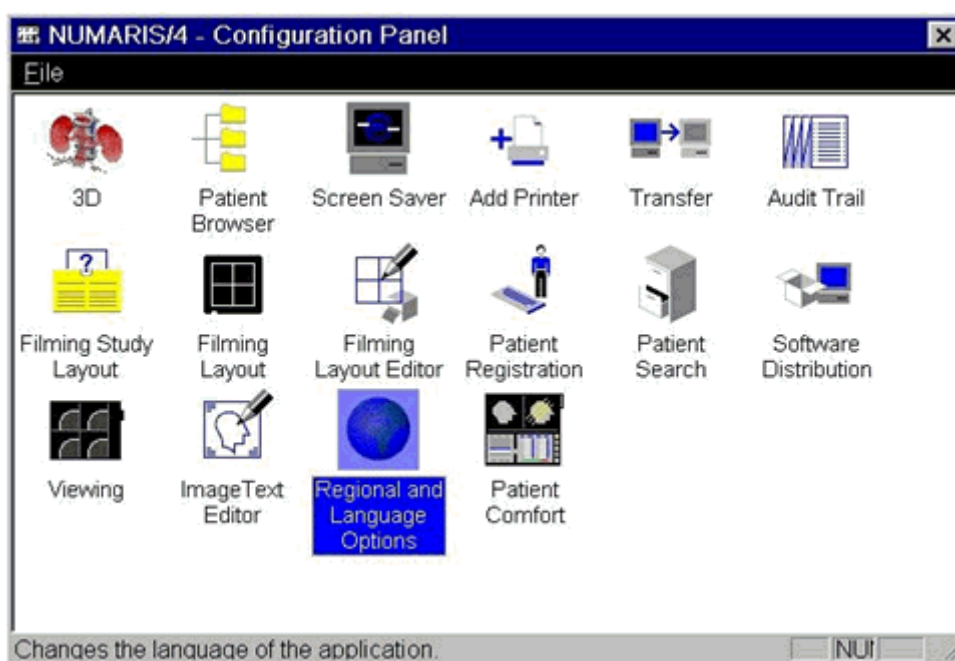

- Select **Japanese (Japan)** in the drop-down menu below "Format"
- Click **Additional settings**

Fig. 208: Regional and Language jp

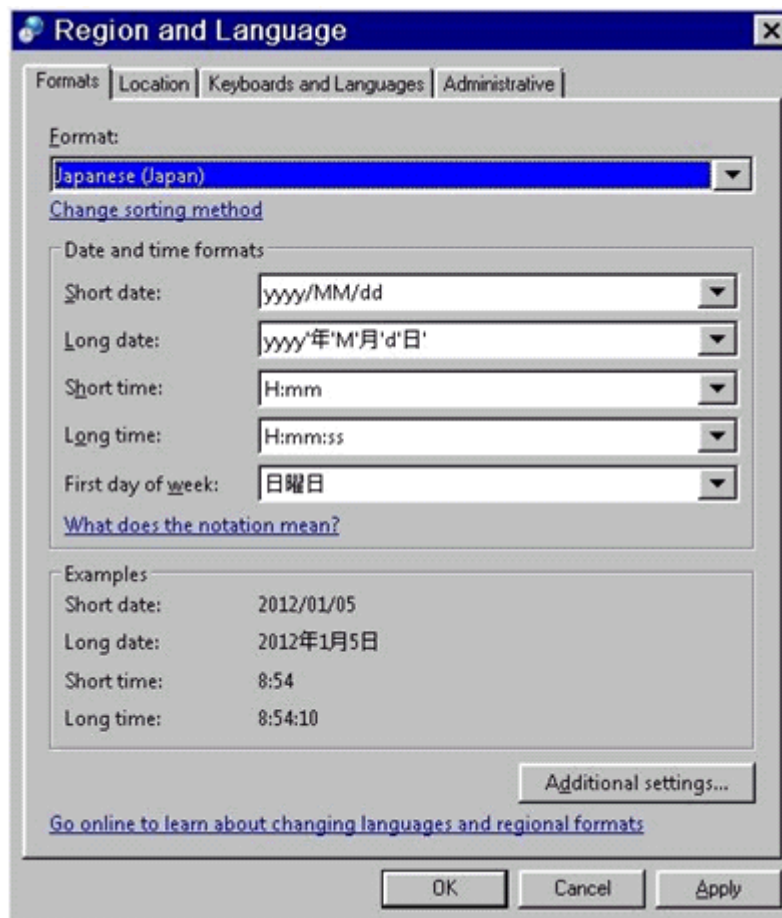


- Select the **Date** tab card
- Select the Calendar drop-down menu

Fig. 209: Calendar type Japan


和暦

Fig. 210: Customize Format - Calendar type

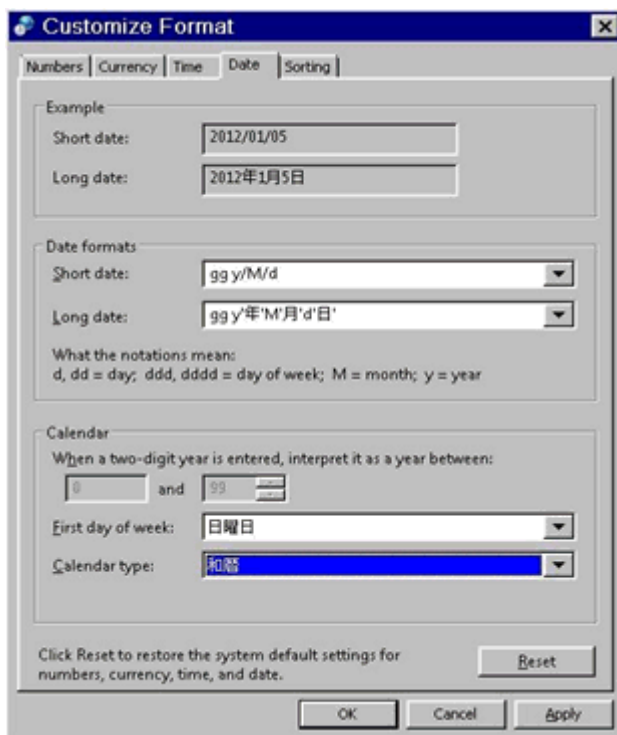

- In the drop-down menu "Current location:" select **Japan**
- Click **OK**

Fig. 211: Region and Language - Location Japan

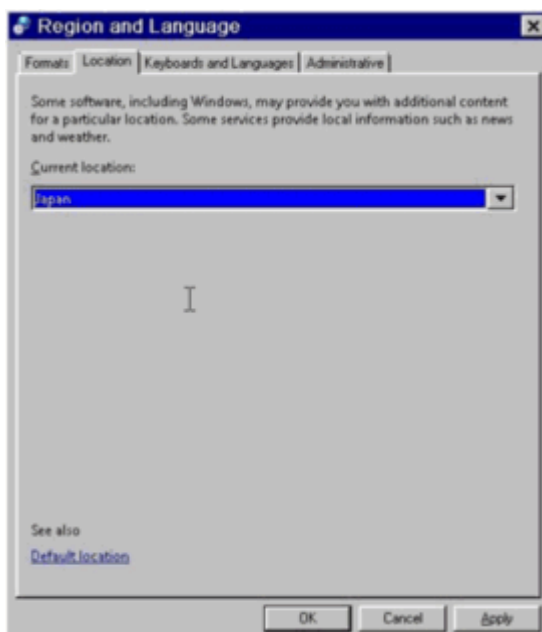


- Select the **Keyboards and Languages** tab card
- In the drop-down menu "Choose a display language:" select:

Fig. 212: Display language Japanese


日本語

Fig. 213: Region and Language - Keyboards and Languages Japan

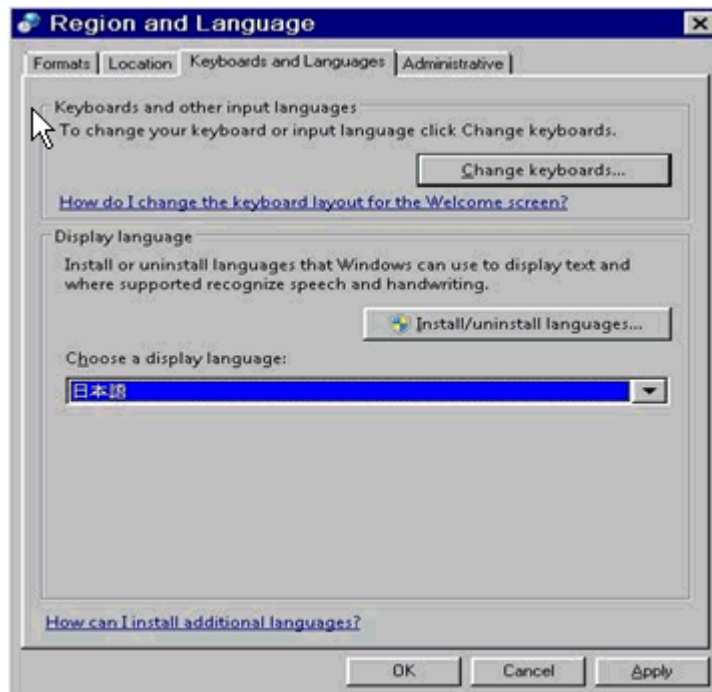


- Select the **General** tab card
- In the drop-down menu "Default input language:" select **Japanese (Japan) - Microsoft IME**:
- Add the Japanese keyboard if not available
- Select **Japanese (Japan) Keyboard Microsoft IME**

Fig. 214: Text Services and Input Languages - General Japan

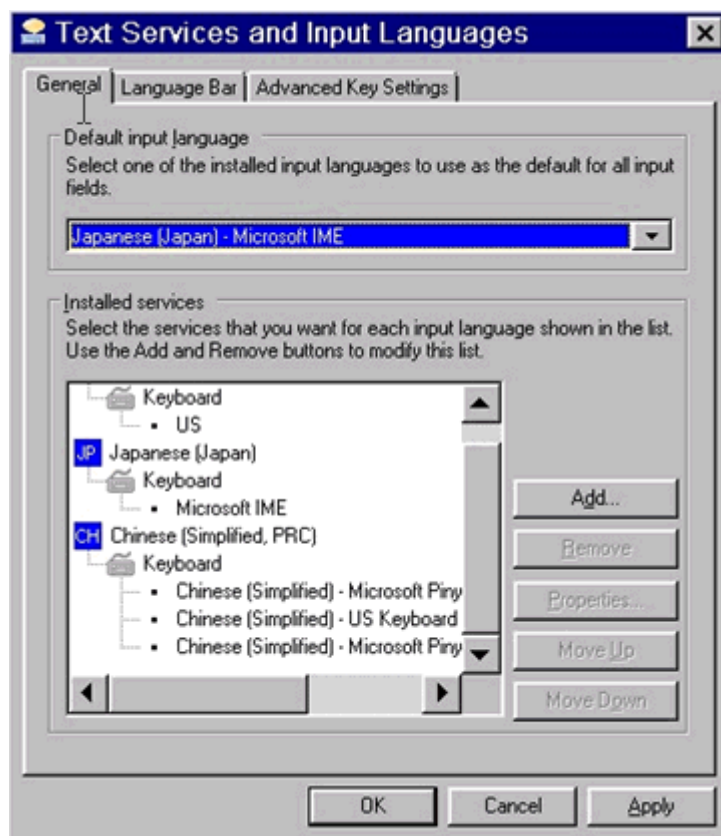

- Set the checkbox **IME** and click **OK**


Fig. 215: Add Input Language Japanese

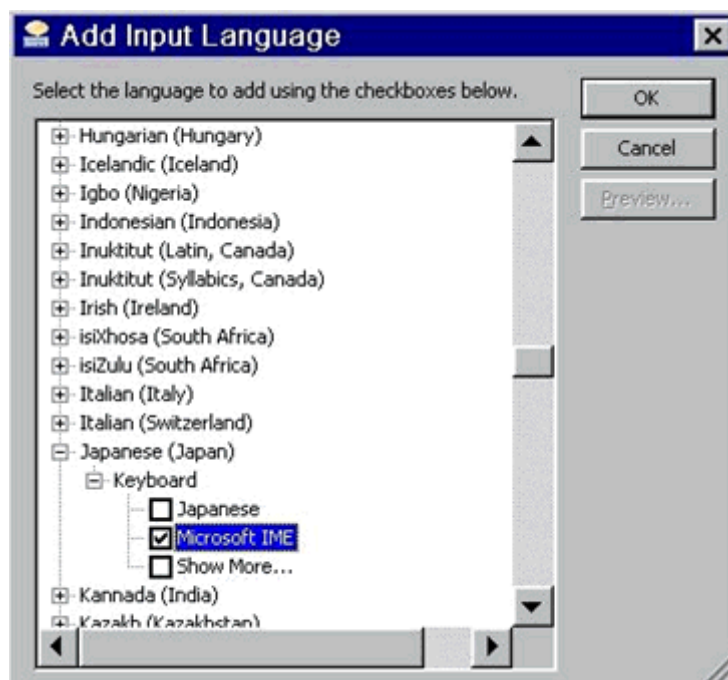

- Select the **Language Bar** tab card
- Select the **Floating On Desktop** radio button


- Set the checkboxes **Show additional Language bar icons in the taskbar** and **Show text labels on the Language bar**
- Click **OK**

Fig. 216: Text Services and Input Languages

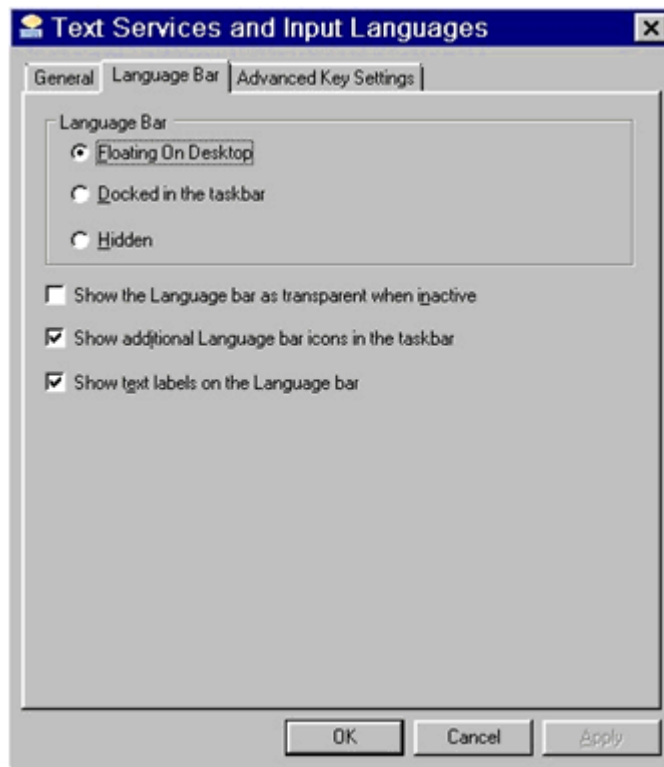

- Select the **Administrative** tab card
- Click **Change system locale**

Fig. 217: Region and Languages - Administrative

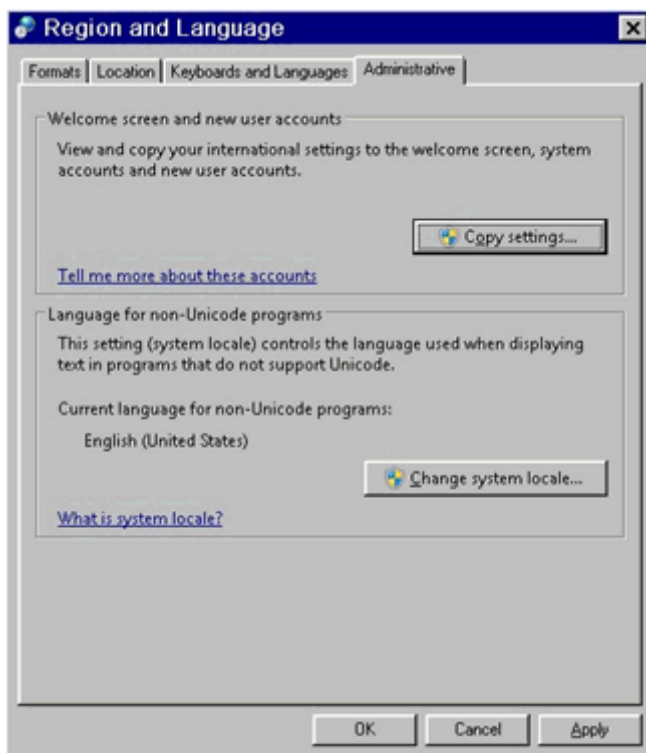


- Click **Apply** in the popup window

Fig. 218: Change Regional Options




- Select Japanese (Japan) in the next popup window
- Click **OK**

Fig. 219: Region and Language Settings - Japanese

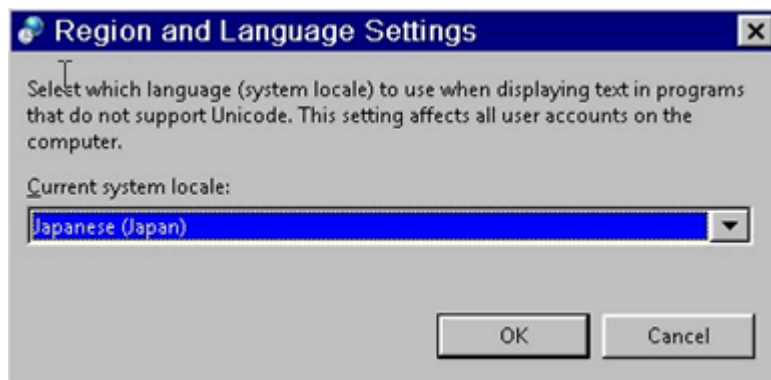


- In the following popup window select **Restart now**

Fig. 220: Change System Locale - Restart now

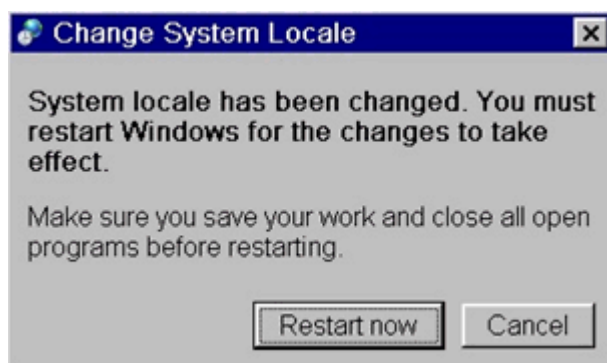


- Confirm this popup with **OK**

Fig. 221: OK - jp

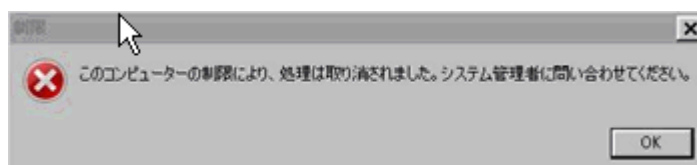

- Confirm the warning with **OK**


Fig. 222: Warning OK - jp

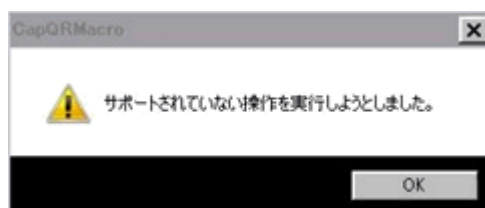

4.1.4.2 Switch to English from Japanese


- Select **Options > Service > Local Service**

Fig. 223: Options - Service - Local Service.jp

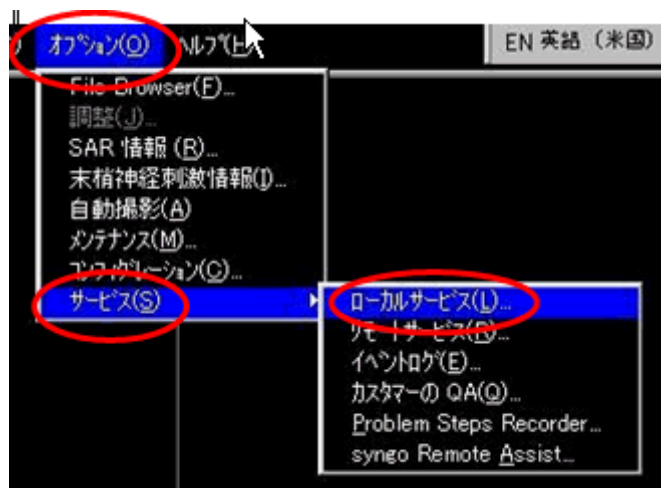

- Select **DICOM Character Set >** in the drop-down menu "DICOM default character set" select **Latin Alphabet No. 1**

Fig. 224: Configuration - Character Set - DICOM default character set (Latin Alphabet No 1).jp

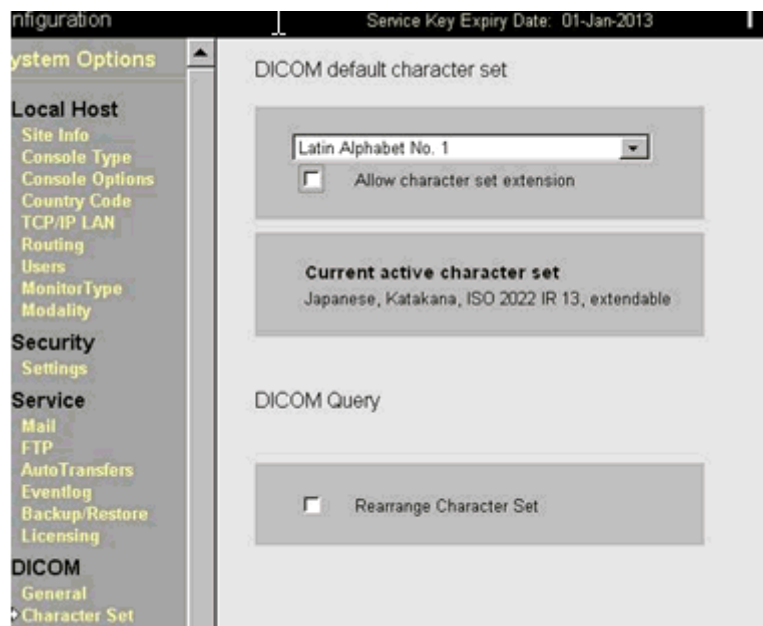


- Confirm "character set data successfully saved" popup with **OK**

Fig. 225: Character set data successfully saved - OK



- Confirm the "application restart popup" with **OK**

Fig. 226: A restart of the application SW is necessary - OK.jp

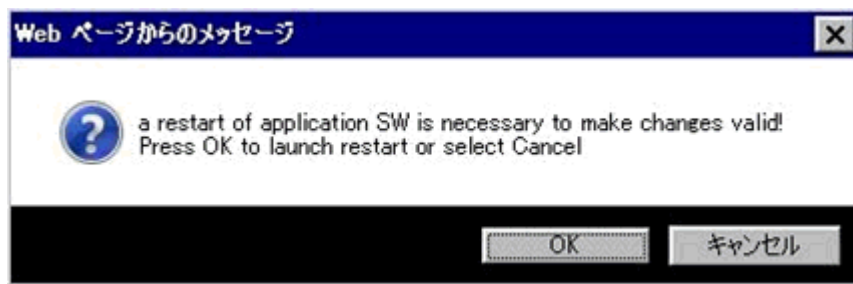


- When the application is up and running again, close the windows "Personalization" and "Local Service"

Fig. 227: Personalization and Local Service - jp




- Click Options > Configuration

Fig. 228: Options Configuration - jp

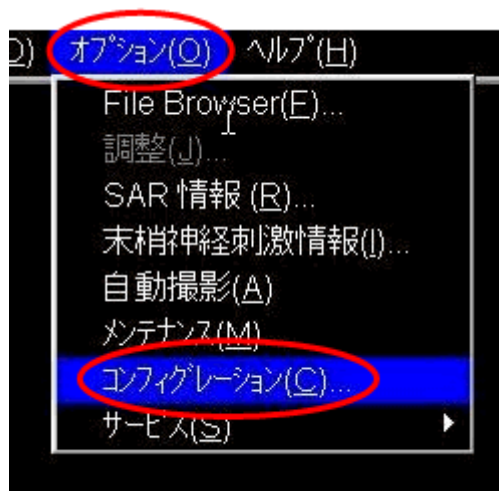

- In the Numaris/4 configuration panel, double-click **Regional and Language Options**


Fig. 229: Regional and Language Options - jp

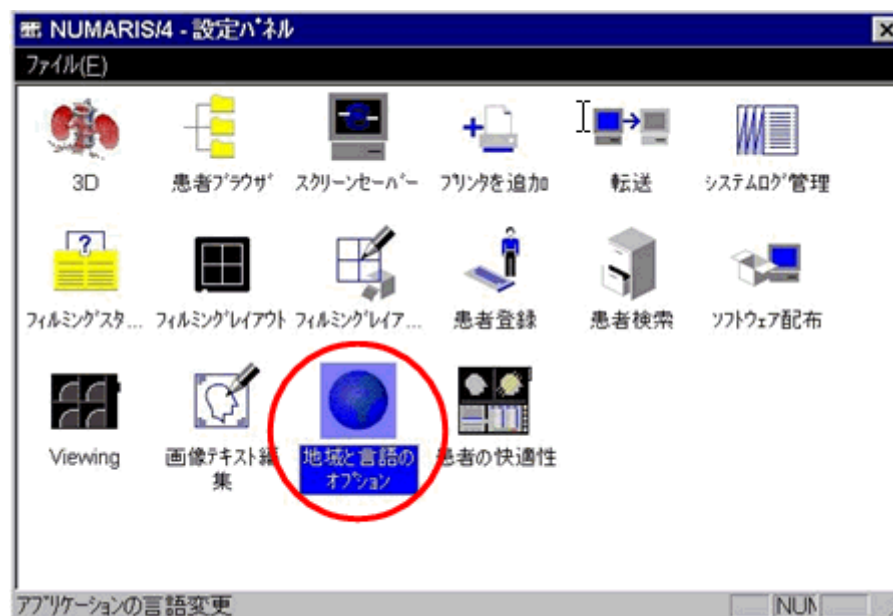

- Set the date time format as follows

Fig. 230: Region and Language - jp


英語 (米国)

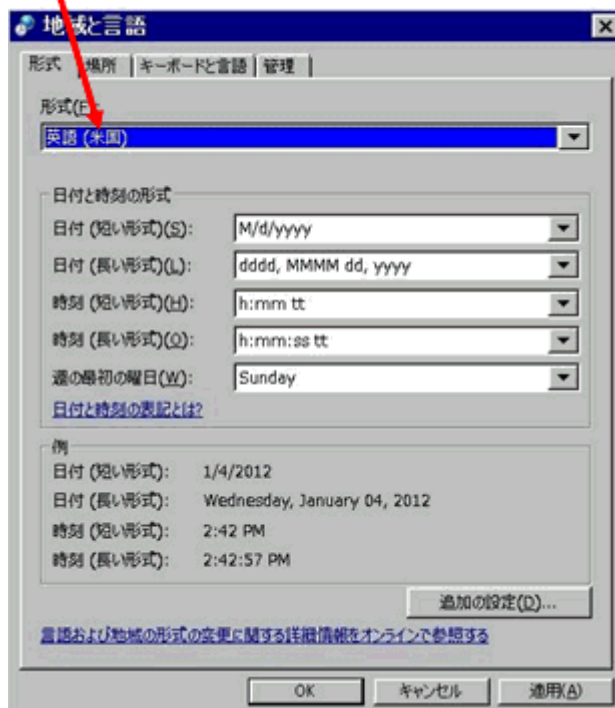

- Select the second tab
- Click OK


Fig. 231: Region and Language - Location jp

米国

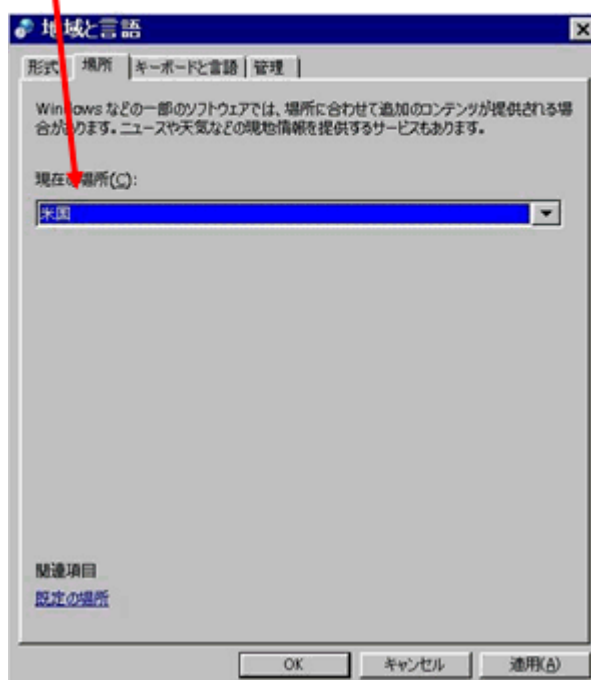

- Select English


Fig. 232: Region and Language - Keyboards and Languages jp

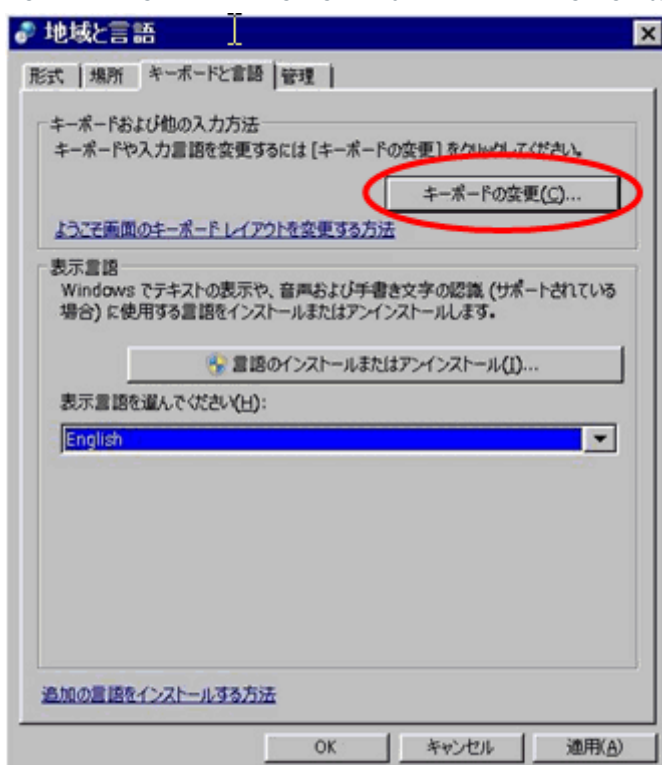

- Click the first tab from the left
- Select the following in the left drop-down menu

Fig. 233: Text Services and Input Languages jp


英語 (米国) - US

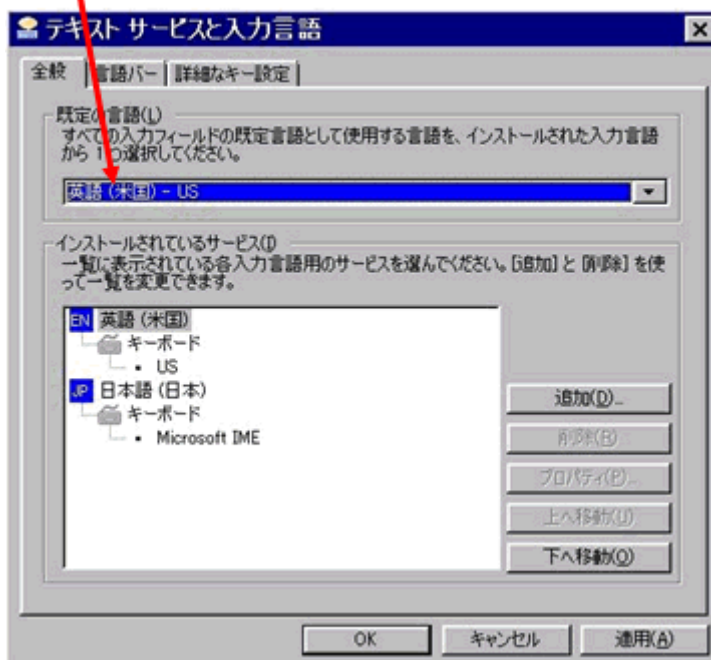

- Select the second tab > select the mid radio button

- Click **OK**

Fig. 234: Text Services and Input Languages jp

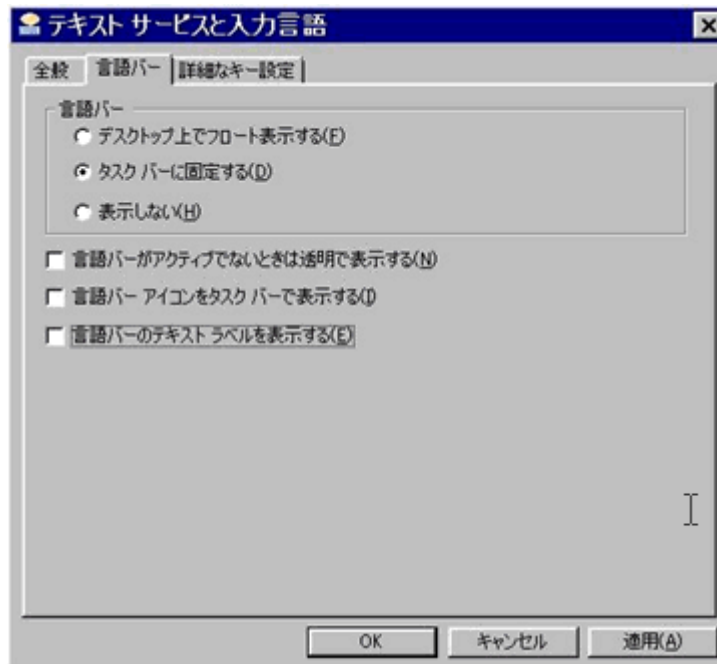

- Select the fourth tab > select the button

Fig. 235: Region and Language - Change system locale jp (zoomed)



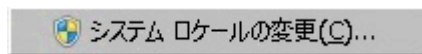


Fig. 236: Region and Language - Change system locale jp

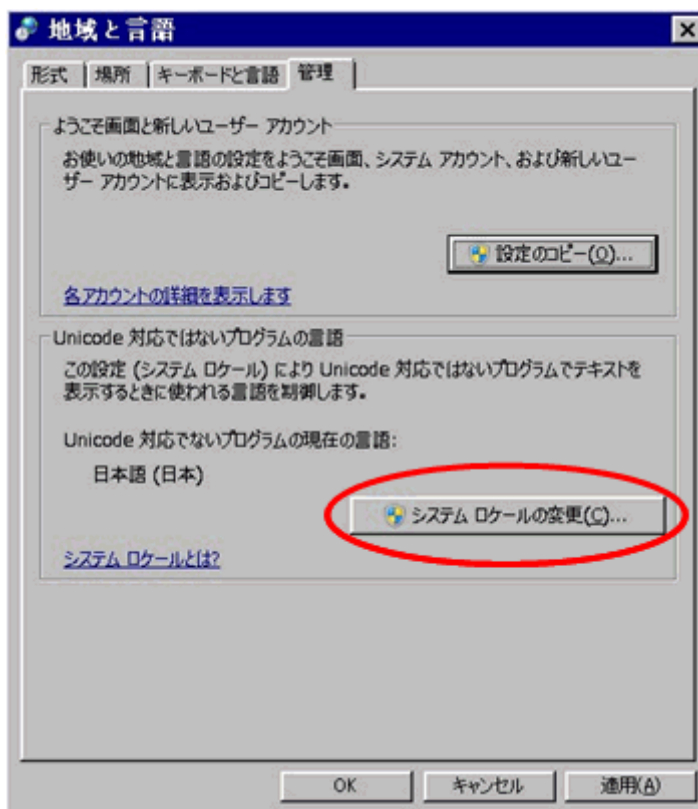


- Click **Apply**.

Fig. 237: Change Regional Options - Apply jp

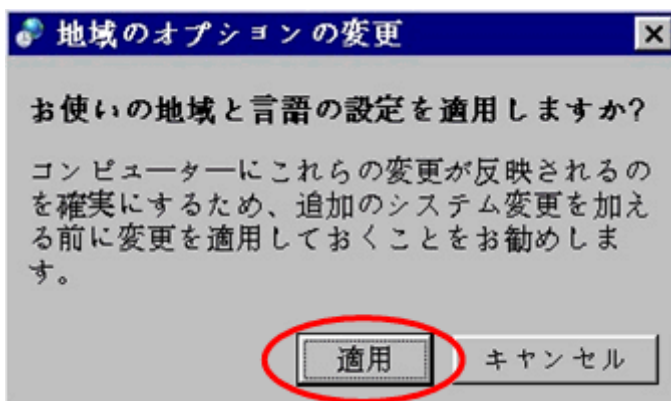


- Select the drop-down menu as follows and click **OK**

Fig. 238: Region and Language Settings - US English jp

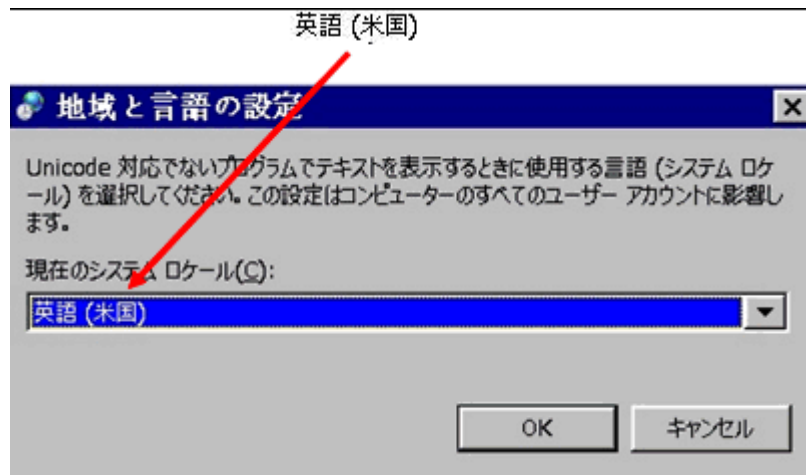

- Proceed as follows:

Fig. 239: Change System Locale - Restart now ch

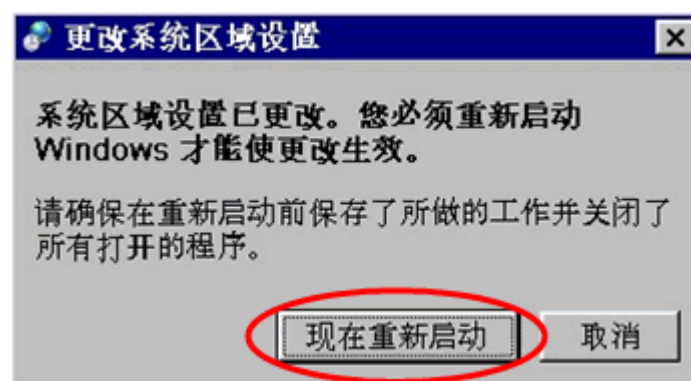

- Confirm the restriction popup with OK


Fig. 240: Confirm Restriction popup with OK



4.1.4.3 Keyboard and Language Settings - Chinese


- In the Service UI, select **Options -> Configuration**
- Select the DICOM **Character Set** > in the drop-down menu "DICOM default character set" select **Chinese**

Fig. 241: Configuration - DICOM default character set

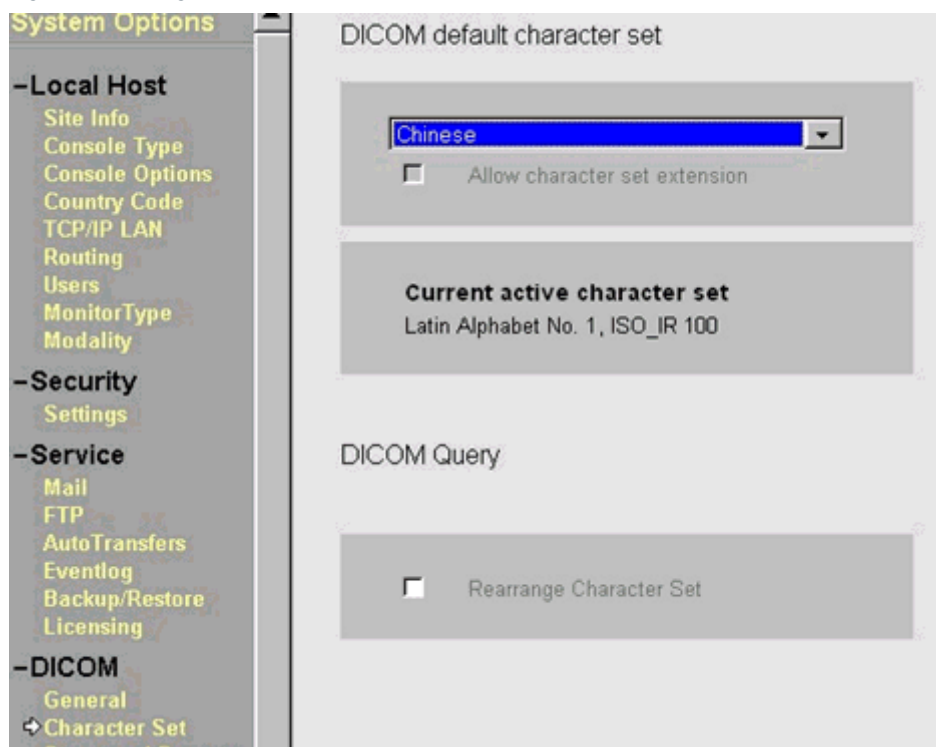

- Click **OK** to confirm the change in character set

Fig. 242: Configuration - character set data successfully saved

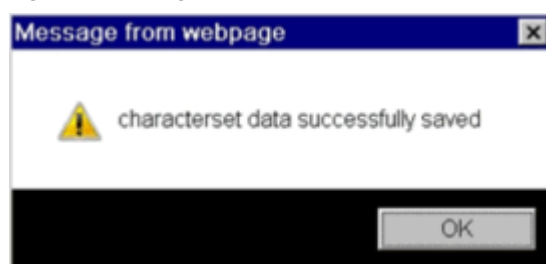


- Click **OK** to restart the application

Fig. 243: Configuration - a restart of the application SW is necessary to make changes valid




- Select **Options > Configuration**
- Select **Regional and Language Options**.

Fig. 244: Configuration panel

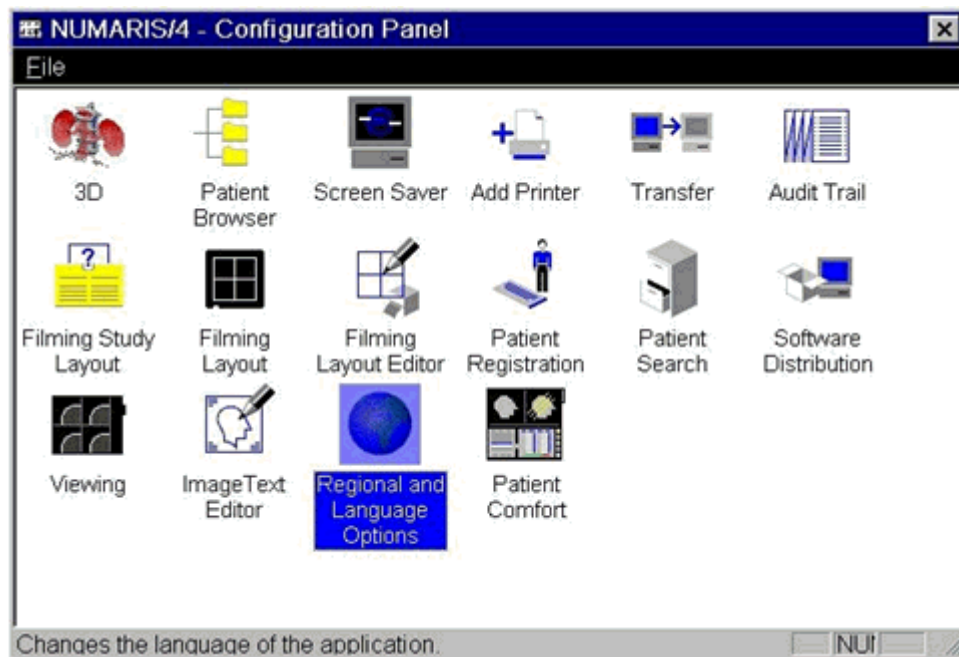


- Select **Chinese (Simplified, PRC)** in the drop-down menu below "Format"

Fig. 245: Regional and Language

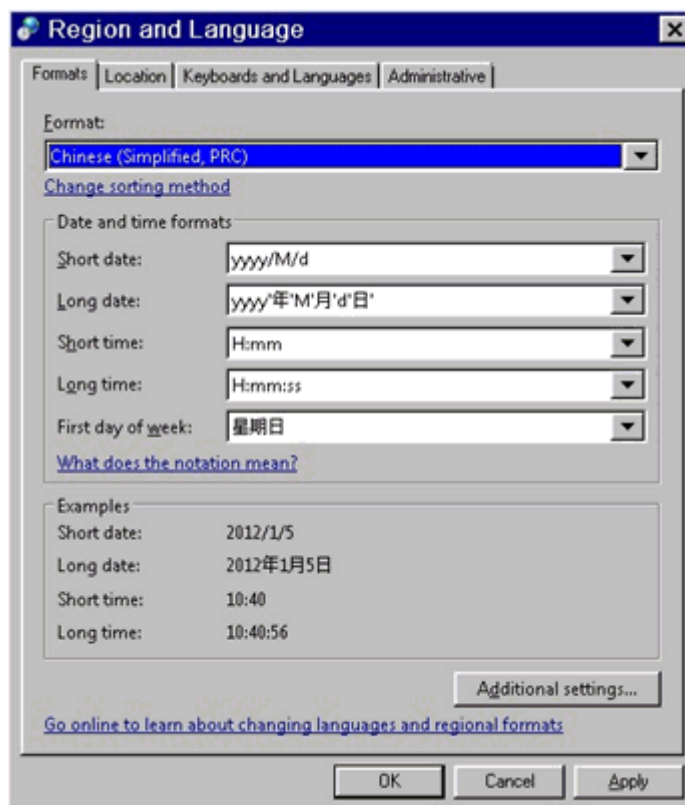

- Select the **Location** tab
- Select **Chinese** in the drop-down menu

Fig. 246: Region and Language - Location China

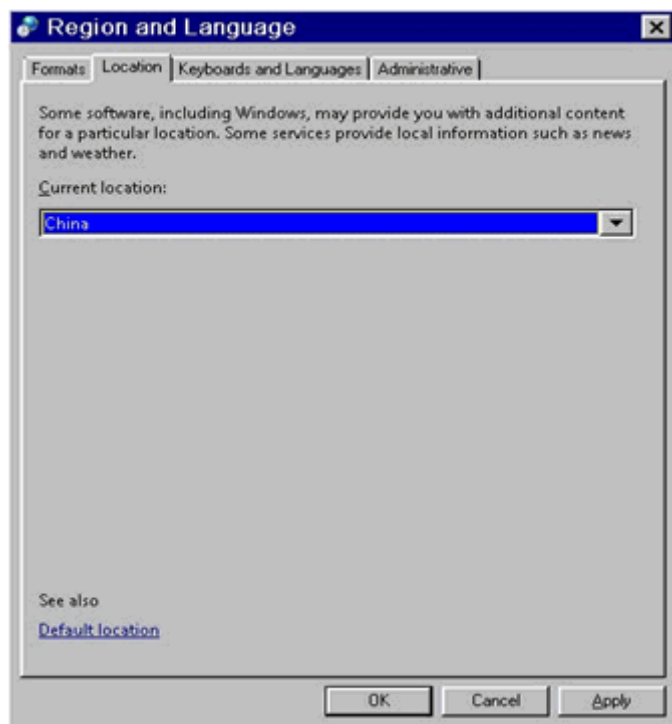


- Select the **Keyboards and Languages** tab
- Select

Fig. 247: Choose display language - Chinese

中文(简体)

- in the drop-down menu **Choose a display language** tab

Fig. 248: Region and Language - Keyboards and Languages

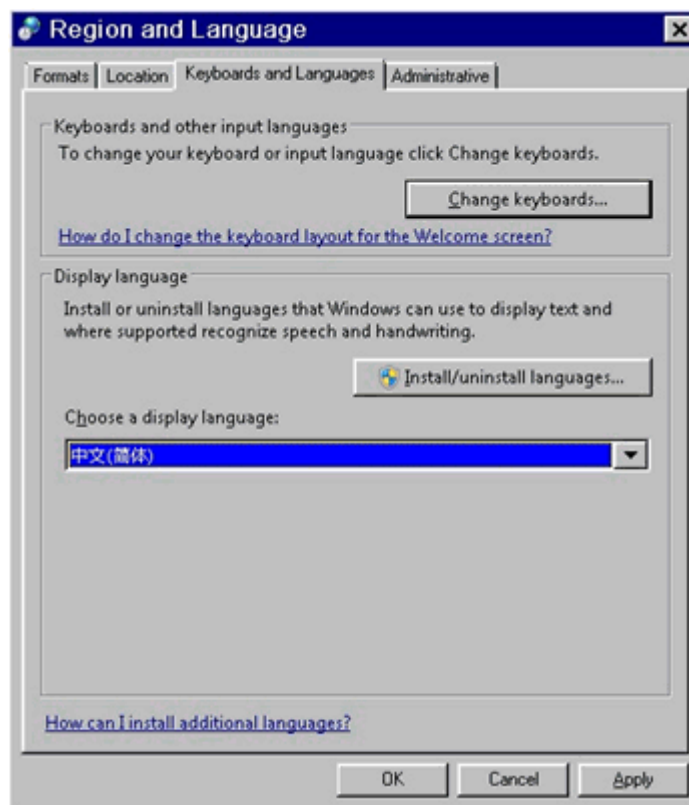


- Select the **General** tab.
- Select the drop-down menu **Chinese (Simplified, PRC) - Chinese (Simplified) - Microsoft Pinyin ABC Input Style**

Fig. 249: Text Services and Input Languages - General




- Set the checkbox “Chinese (Simplified, PRC) - Chinese (Simplified) - Microsoft Pinyin ABC Input Style”

Fig. 250: Add Input Language

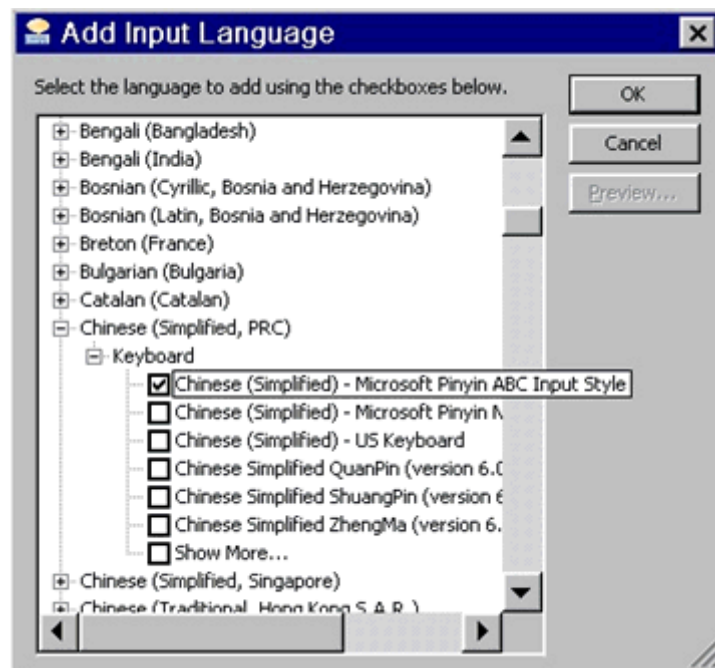


- Select the **Language Bar** tab card
- Select the **Floating On Desktop** radio button
- Set the checkbox **Show additional language bar icons in the task bar**
- Set the checkbox **Show text labels on the language bar**

Fig. 251: Text Services and Input Languages

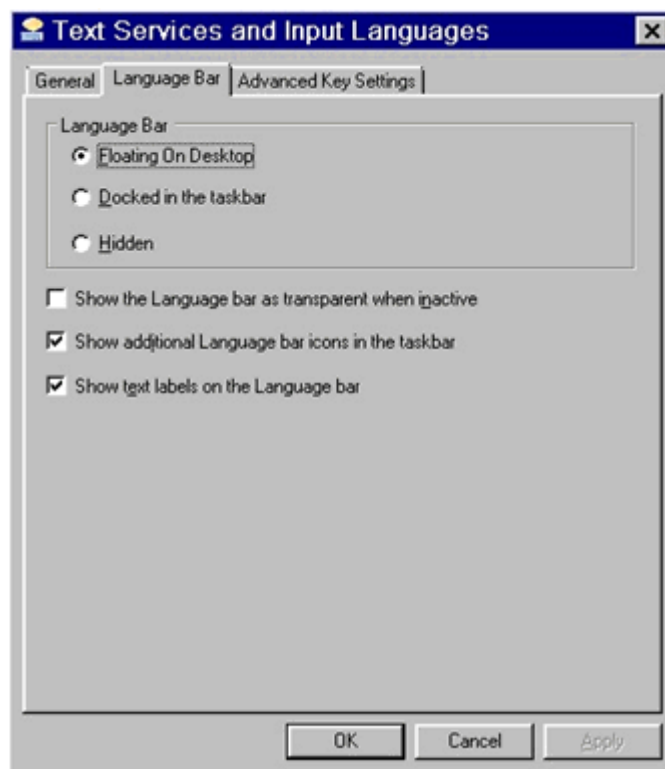

- Select the **Administrative** tab
- Click **Change system locale**

Fig. 252: Region and Languages - Administrative

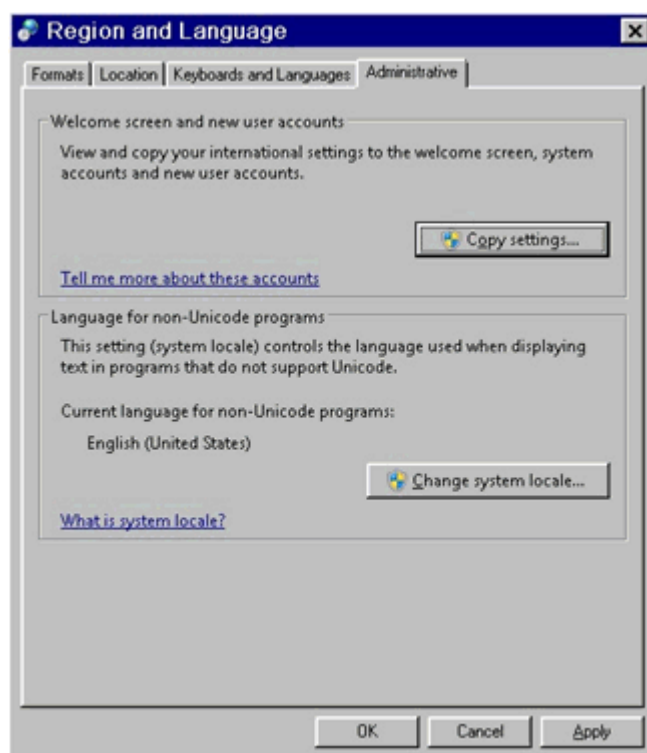

- Click **Apply** to confirm the region and language changes

Fig. 253: Change Regional Options




- Select **Chinese (Simplified, PRC)** in the drop-down menu "Current system locale"

Fig. 254: Region and Language Settings - Chinese




- Click **Restart now**

Fig. 255: Change System Locale - Restart now

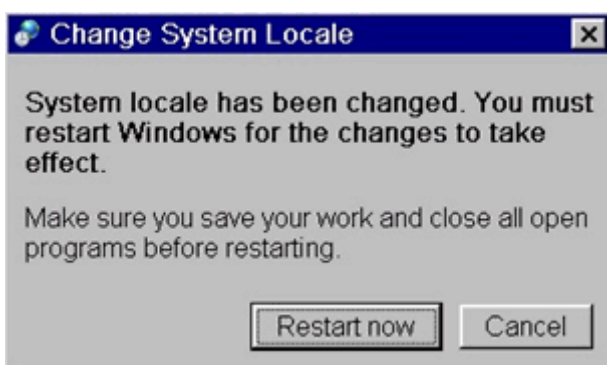

- Log in as meduser, wait until the application is up and running
- Confirm the following popups with **OK**

Fig. 256: OK - ch

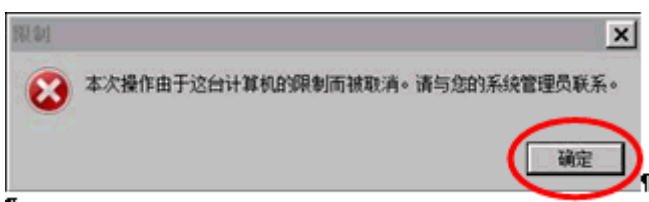

- Confirm the warning with **OK**


Fig. 257: Warning OK - ch



4.1.4.4 Switch to English from Chinese


- In the Service UI, select **Options** -> **Configuration**
- Select the DICOM **Character Set** > in the drop-down menu "DICOM default character set" select **Chinese**
- Open the Local Service service GUI

Fig. 258: Options - Service - Local Service ch

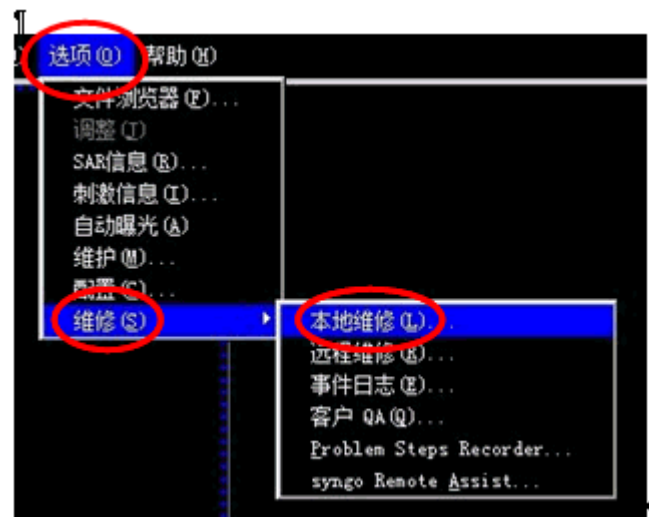


- Set the DICOM Character Set to **Latin Alphabet No. 1**

Fig. 259: Configuration - Character Set - DICOM default character set (Latin Alphabet No 1)

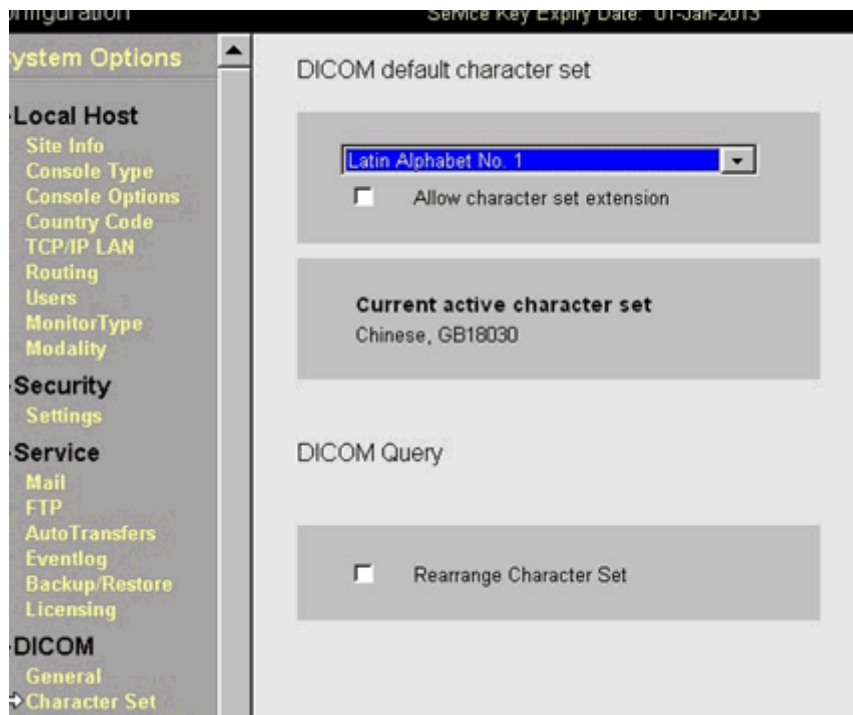

- Confirm the "character set data successfully saved" popup with **OK**

Fig. 260: Character set data successfully saved - OK ch

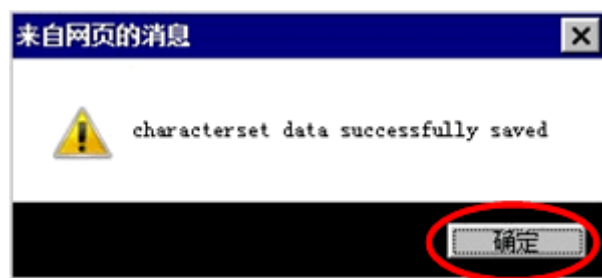

- Confirm the application restart popup with **OK**

Fig. 261: A restart of the application SW is necessary - OK ch




- When the application is up and running again, close the windows "Personalization" and "Local Service"

Fig. 262: Personalization and Local Service - ch



- Click **Options > Configuration**

Fig. 263: Options Configuration - ch

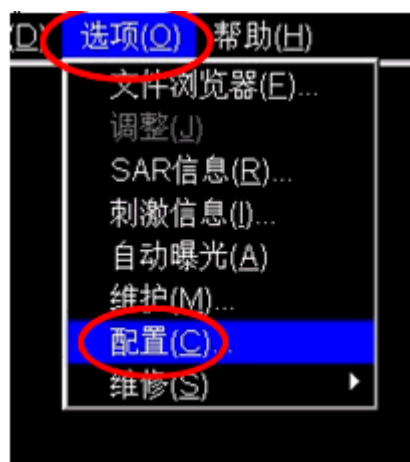


- In the Numaris/4 configuration panel, double-click **Regional and Language Options**

Fig. 264: Regional and Language Options - ch




- Set the date time format as follows

Fig. 265: Region and Language - ch

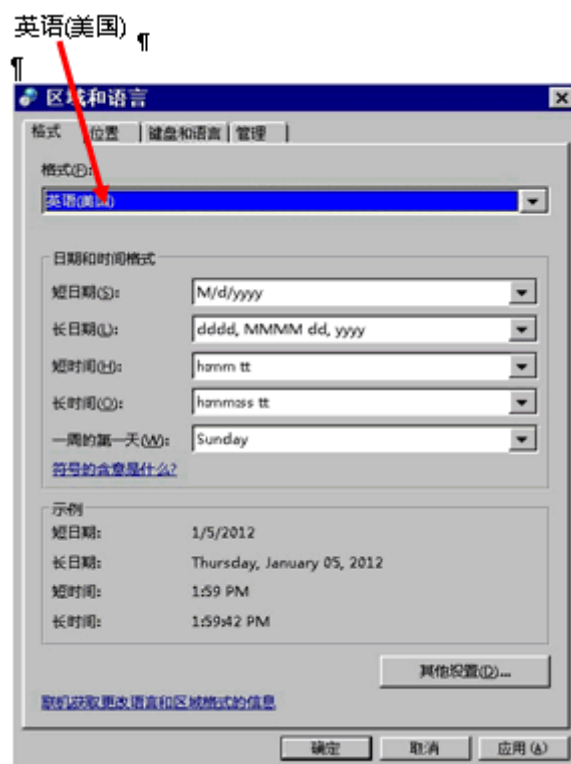


- Select the second tab
- Click OK.

Fig. 266: Region and Language - Location ch

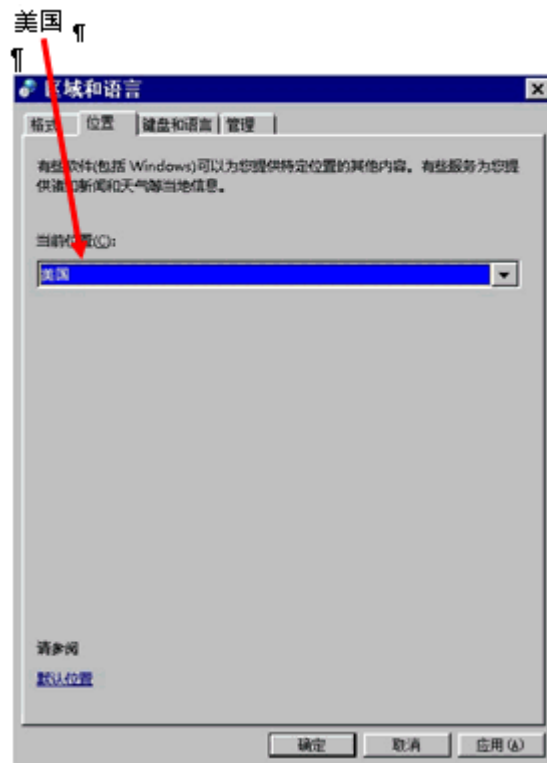


- Select **English**

Fig. 267: Region and Language - Keyboards and Languages ch

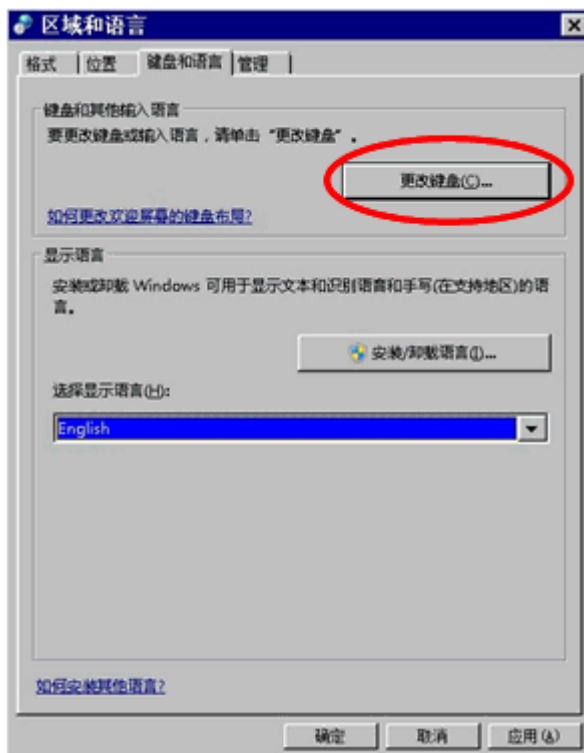


- Click the first tab from the left
- Select the following in the left drop-down menu

Fig. 268: Text Services and Input Languages ch

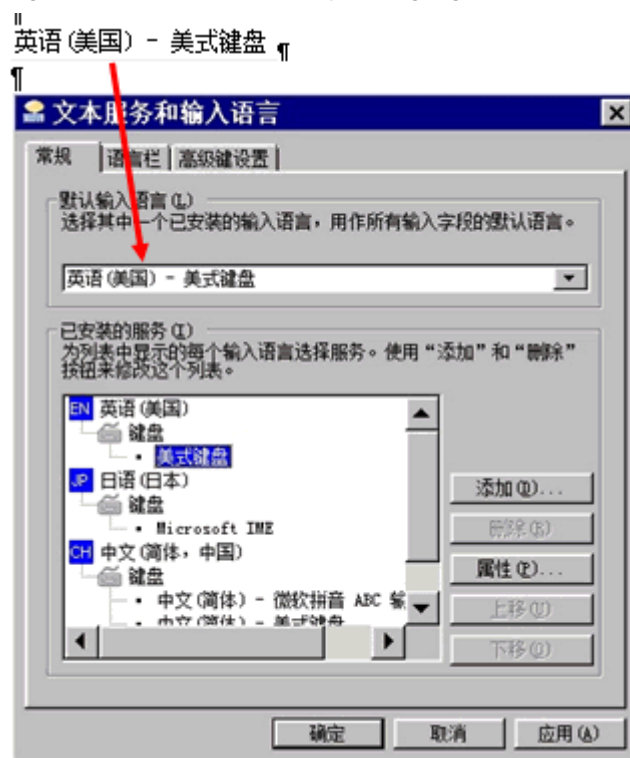

- Select the second tab > select the mid radio button
- Click OK

Fig. 269: Text Services and Input Languages ch

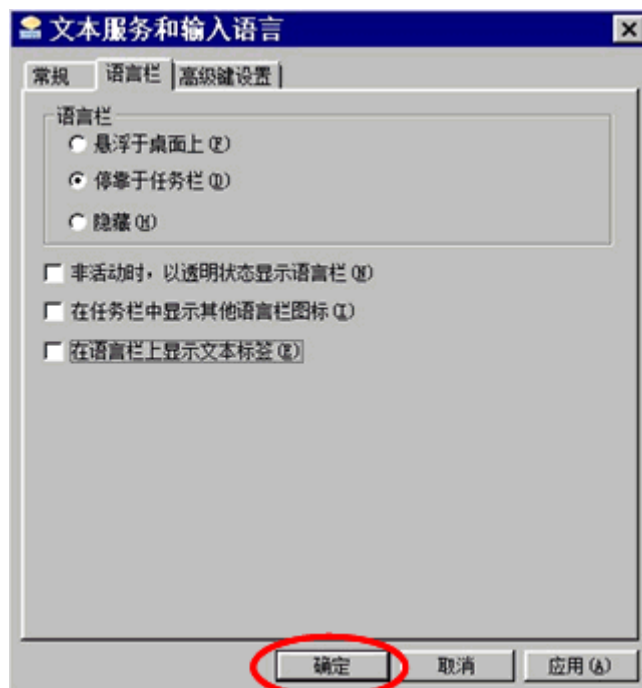


- Select the fourth tab > select the button

Fig. 270: Region and Language - Change system locale ch (zoomed)




Fig. 271: Region and Language - Change system locale ch

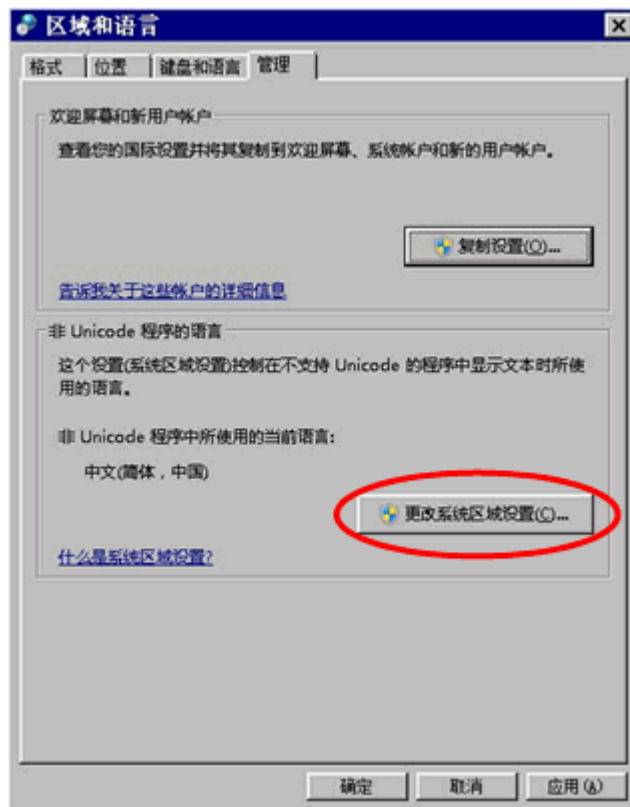


- Click **Apply**.

Fig. 272: Change Regional Options - Apply ch

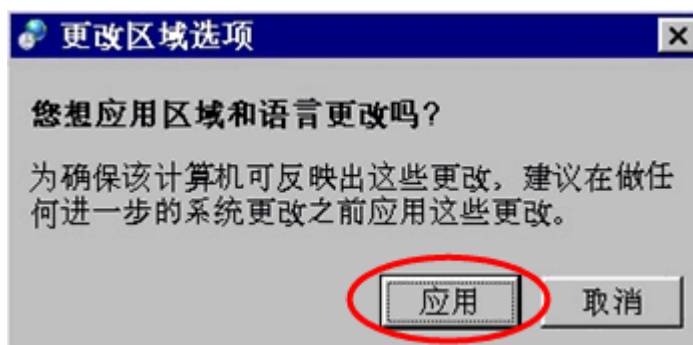


- Select the drop-down menu as follows and click **OK**

Fig. 273: Region and Language Settings - US English ch

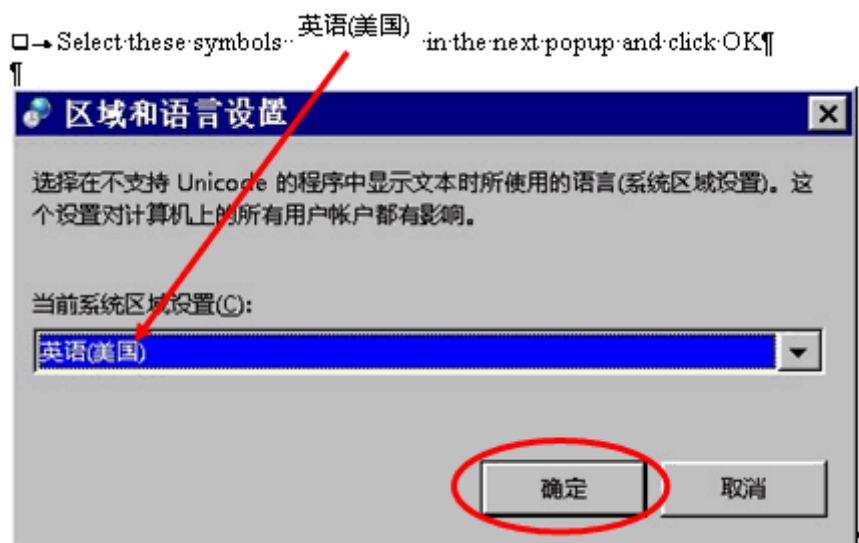

- Proceed as follows:

Fig. 274: Change System Locale - Restart now ch

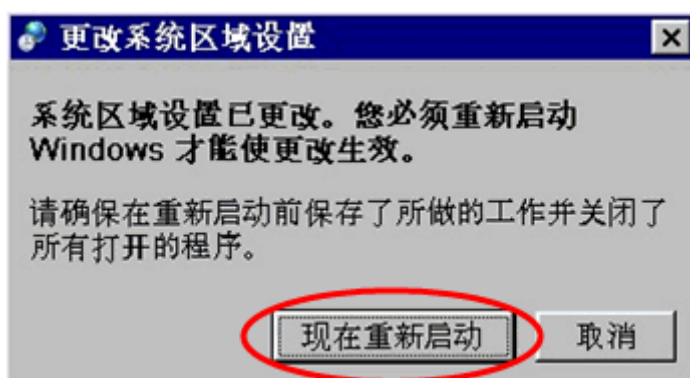

- Confirm the restriction popup with OK


Fig. 275: Confirm Restriction popup with OK




4.1.5 Deactivation of OLM

The deactivation of the OLM (optical length meter) will improve system reliability. To ensure the correct travel distance of the tabletop and prevent system error messages caused by the deactivated OLM, it is very important that the compatible system software version is installed immediately before performing the following steps. The compatible system software is needed to enter the table correction factor; non-compatible system software does not show the correction factor in the CMD tool.

4.1.5.1 Preparation

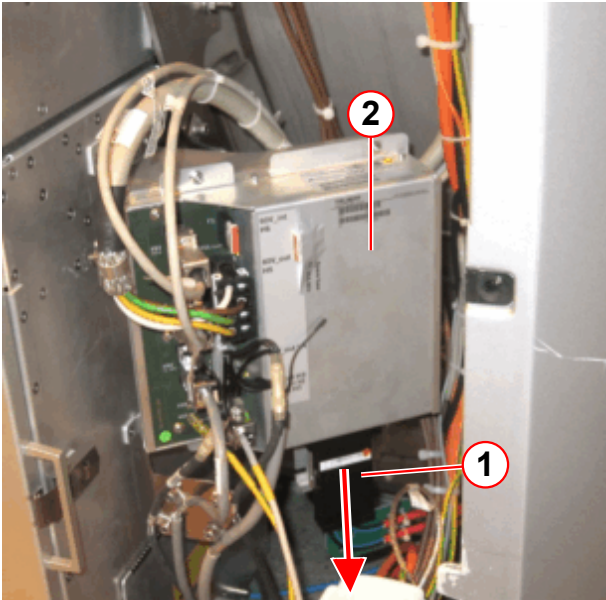ACTION	MRAWP
Check the currently installed system software version, the compatible system software version must be installed N4_VE11A_LATEST_20140830_P12 (SP01) before continuing with these steps.	☐
Move the patient table into the Home position.	☐
Shut down and switch off the system.	☐
Remove the service access cover on the left side of the magnet to access the OLM plug (XB5).	☐
The OLM plug is connected to the PAM at the right side of the RFCEL. Loosen the two securing screws for OLM plug XB5 and pull off the partially magnetic plug (1). <i>Fig. 276: Pulling off OLM plug XB5.</i>  (1) OLM plug XB5 (2) PAM:	☐

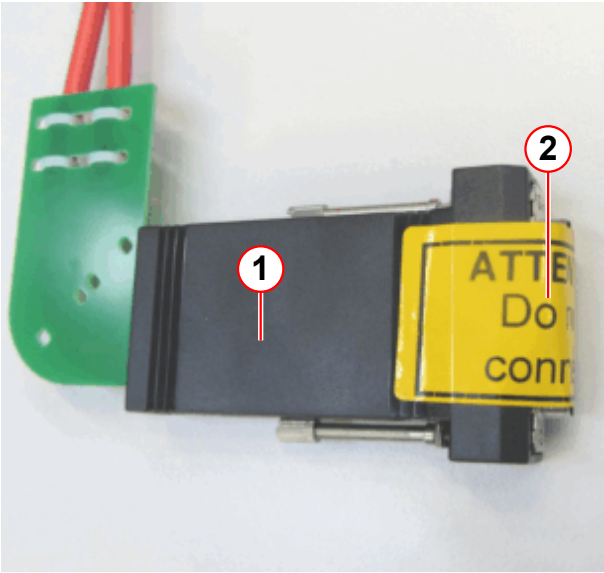
ACTION	MRAWP
Affix the delivered self-adhesive label "Attention! Do not connect!" to OLM plug XB5 (2). <i>Fig. 277: Labeling OLM plug XB5.</i>  (1) OLM plug XB5 (2) "Attention! Do not connect!" label	☐
Affix the OLM plug to the cable harness next to the PAM with a cable tie. <i>Fig. 278: Setting the location</i> 	☐
Switch on the system.	☐


NOTICE



The following work steps must be done for each fixed and mobile patient table.


- » If the system has several mobile tables (e.g. second mobile table), repeat the work steps for each mobile patient table.

4.1.5.2 Check of the tabletop travel distance

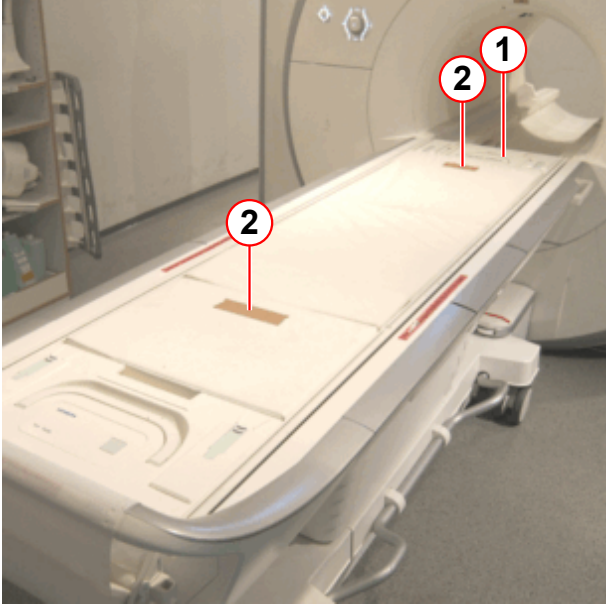ACTION	MRAWP
■ Remove everything except the cushions from the tabletop.	☐
■ Affix pieces of the delivered tape approximately centered on the tabletop cushions (2), with a distance of approx. 30 cm and 170 cm measured from the head end of the tabletop (1). <i>Fig. 279: Preparing the tabletop</i>  (1) Tabletop head end (2) Tape at a distance of approx. 30 cm and 170 cm from the head end.	☐

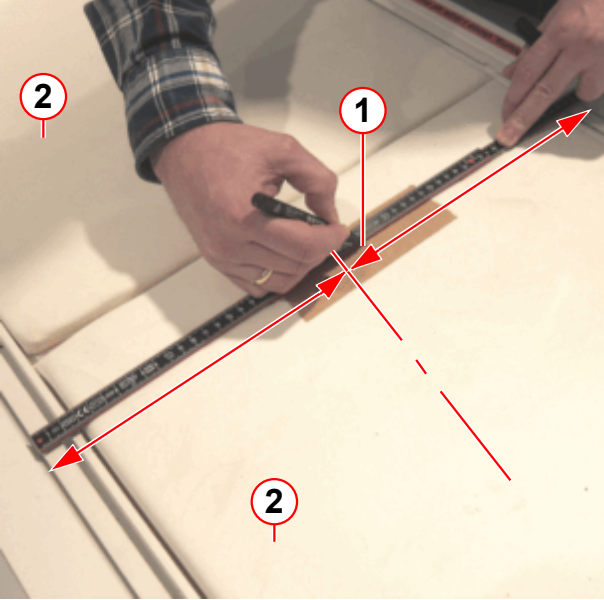
ACTION	MRAWP
Measure the width of the tabletop and mark the center on the two affixed pieces of tape. <i>Fig. 280: Marking the tabletop center</i>  (1) Affixed tape positioned approx. 30 cm and 170 cm from the head end of the tabletop (2) Cushions on tabletop	<input data-bbox="1353 275 1393 309" type="checkbox"/>

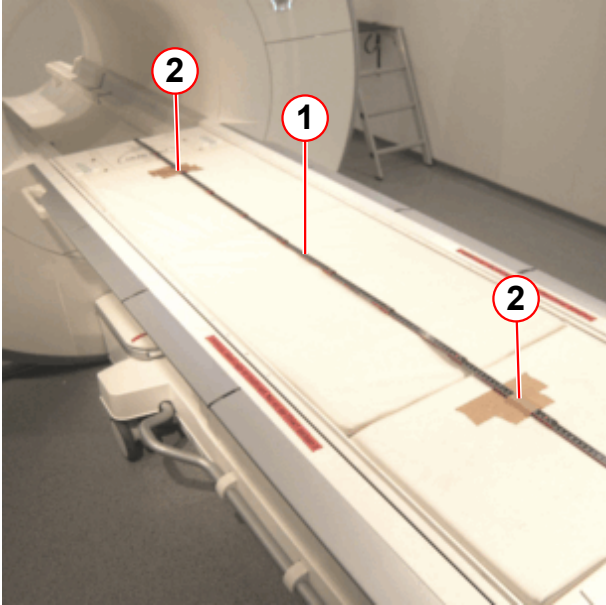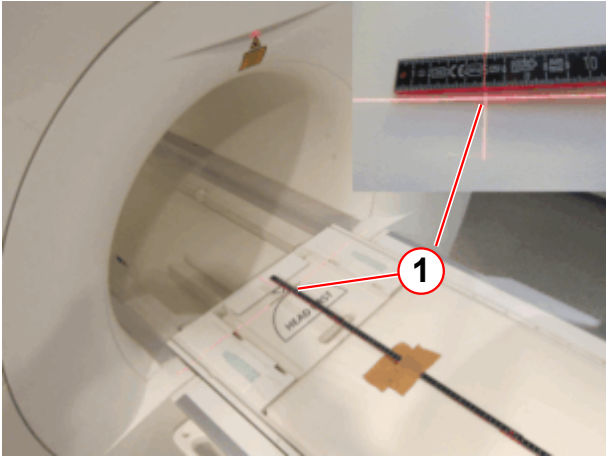ACTION	MRAWP
Position the zero point of the folding ruler scale at the head end of the tabletop. Align the 2 m long folding ruler (1) with the center marks and affix the ruler with two pieces of tape (2). Ensure that the folding ruler is completely straight and has no tension. Use only the EU standard accuracy class 2 folding ruler (p. n. 10759380) delivered with the update kit. Accuracy class 2 can be identified by the sign II inside an oval circle on the ruler. <i>Fig. 281: Positioning the folding ruler</i>  (1) Folding ruler (2) Tape and center marks	☐
Switch on the patient table light marker and position the tabletop at exactly 50 mm from the folding ruler. <i>Fig. 282: Positioning at 50 mm</i>  (1) Light marker at 50 mm position of folding ruler	☐
Switch off the light marker.	☐

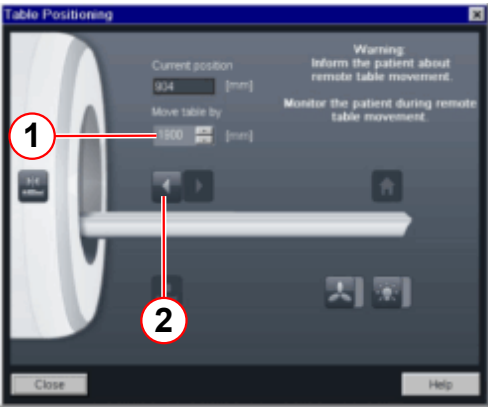ACTION	MRAWP
Open the "Table Positioning" window on the MR console. <i>Fig. 283: Move table by 1900 mm.</i>  (1) Move table by (entered value of 1900 mm). (2) Start button for table movement towards the magnet.	☐
Enter 1900 (1) for the "Move table by" setting.	☐
Start table movement into the magnet by clicking the arrow toward the magnet (2).	☐
Switch on the light marker and read the distance value at the folding ruler. The value must be 1950 +/- 2 mm. Record the distance value: Distance = _____	☐
Continue with the next work steps only if the travel distance is out of specification. If the travel distance is within specification skip the next working steps and continue with Labeling(→ Labeling the horizontal drive gear box assembly / Page 226).	☐

4.1.5.3 Open the CMD tool window

ACTION	MRAWP
Select Options, ⇒ Service, ⇒ Local Service.	☐
Enter the valid service key (service level 7) and click OK.	☐
Select Utilities in the service platform.	☐

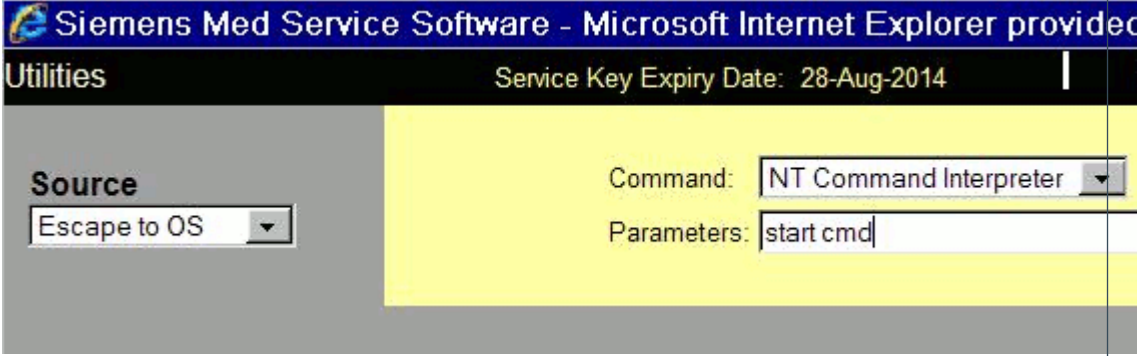ACTION	MRAWP
■ Start the CMD tool window from the local service UI under Utilities: - Under Source, select <code>Escape to OS</code>, - Under Command, select <code>NT Command Interpreter</code> - Under Parameters, type in the following: <code>Start CMD</code> and press Go to execute the program. <i>Fig. 284: Start CMD tool window from local service</i> 	☐
■ A CMD tool window is displayed	☐



4.1.5.4 Calculation of the new table correction factor (only needed if the measured travel distance is out of specification)

ACTION	MRAWP
Select the cmd tool windows, enter the case-sensitive command MrVerifyPTAB ↵, and record the actual correction factor (default = 100150): CorrFactor (actual) = _____	☐
Calculate the new table correction factor with the following formula: CorrFactor (new) = CorrFactor (actual) - Y Y = 50 x travel distance correction in mm. For example, if the travel distance is 3 mm too short: if the actual correction factor is 100150 and the travel distance is 3 mm too short, the calculation would be: CorrFactor (new) = 100150 - (50 x 3) ⇒ The new correction factor in this case would be 100150 - 150 = 100000. Travel distance 3 mm too long: if actual correction factor is 100150 and the travel distance is 3 mm too long, the calculation would be: CorrFactor (new) = 100150 - [50 x (-3)] ⇒ The new correction factor in this case would be 100150 + 150 = 100300. Record the calculated CorrFactor (new) = _____	☐

4.1.5.5 Entering the new table correction factor


ACTION	MRAWP
Press the table stop button on the intercom operating unit (this will switch off the table controller, which is necessary for the correction factor change).	☐
Select the cmd tool window and enter the calculated CorrFactor (new) (in this example 100300) with the following command (case-sensitive command with blanks): MrVerifyPTAB -cF 100300 -correction ↵	☐
Enter the command MrVerifyPTAB ↵ and recheck the new correction factor value.	☐
Reset the table stop error on the patient table operating panel.	☐

4.1.5.6 Recheck of the tabletop travel distance

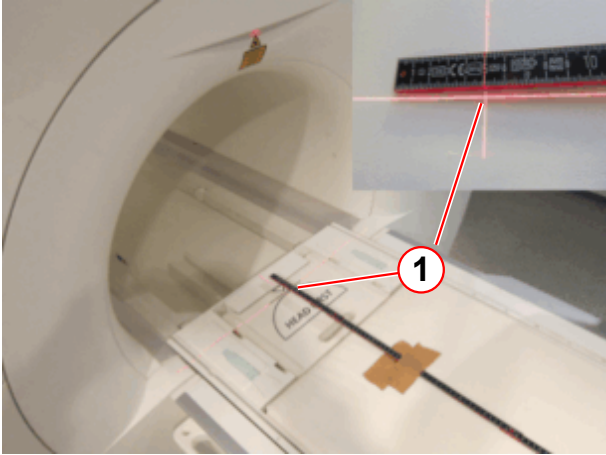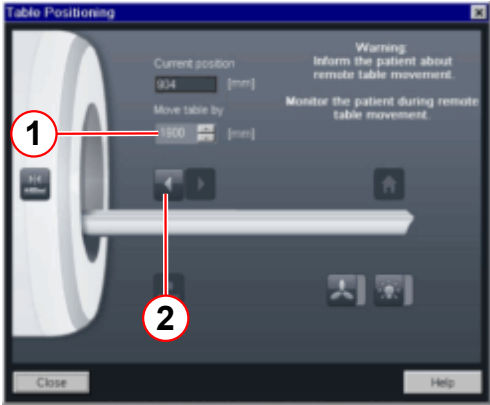ACTION	MRAWP
Switch on the patient table light marker and position the tabletop at exactly 50 mm from the folding ruler. <i>Fig. 285: Positioning at 50 mm</i>  (1) Light marker at 50 mm position of folding ruler	☐
Switch off the light marker.	☐
Open the "Table Positioning" window on the MR console. <i>Fig. 286: Move table by 1900 mm.</i>  (1) Move table by (entered value of 1900 mm). (2) Start button for table movement towards the magnet.	☐
Enter 1900 (1) for the "Move table by" setting.	☐
Start table movement into the magnet by clicking button (2).	☐

ACTION	MRAWP
Switch on the light marker and read the distance value at the folding ruler. The value must be 1950 +/- 2 mm. Record the distance value: Distance = _____	☐
If the travel distance is still out of specification, calculate and enter the new correction factor as described. Repeat this procedure until the travel distance is within specification. If the travel distance is within specification continue with (→ Labeling the horizontal drive gear box assembly / Page 226).	☐

4.1.5.7 Labeling the horizontal drive gear box assembly

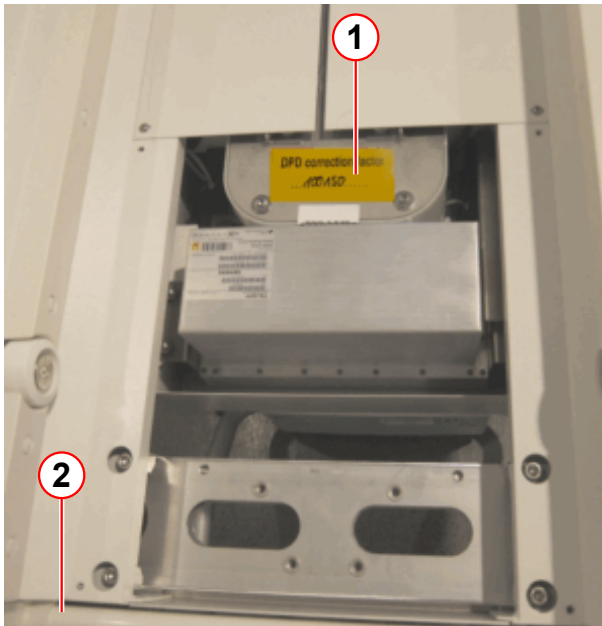ACTION	MRAWP
With a permanent marker write the table correction factor (new) on the delivered "DPD correction factor" self-adhesive label. For tables where the correction factor is unchanged, write the default value 100150 on the label.	☐
Drive the tabletop about 1 m into the magnet.	☐
Remove the cover at the foot end of the patient table (2) and affix the "DPD correction factor" label to the horizontal drive gear box (1).	☐
<i>Fig. 287: Labeling the horizontal drive gear box assembly</i>  (1) New correction factor label at horizontal drive gear box assembly (2) Table foot end	
Affix the "DPD correction factor" label, marked with the corresponding table correction factor, to each patient table used with the system (second mobile table).	☐
Install the foot end table cover.	☐
Close the table positioning window and the CMD tool window.	☐

ACTION	MRAWP
■ Remove the folding ruler and the tape from the tabletop.	☐
■ Close the service access cover.	☐

4.1.6 New Password Handling in syngo



There are no longer fixed or default passwords within syngo. Only random or individually defined passwords are valid. OS Hardening is enabled with Numaris VE11C. With OS Hardening enabled, a specific Administrator account has to be defined. The default "SYADMIN" is suggested. This account is created during the installation in the factory. Autologin for the meduser is also enabled during the factory installation.

4.1.6.1 Disable Autologin

ACTION	MRAWP/MRWP
■ Select Options ⇒ Service ⇒ Local Service .	☐
■ Enter the valid service key (service level 4 or 7) and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select Local Host Users	☐
■ Remove the Enable Autologin checkbox mark to disable the Autologin function.	☐

Fig. 288: SESO users



ACTION	MRAWP/MRWP
■ The password for both users (meduser and syadmin) must now be set. The password must meet the following conditions: - Minimum length of 14 characters - It should have 2 uppercase letters and 2 lowercase letters - It should have 2 numbers and 2 special characters.	☐
■ Click Save .	☐
■ Click Finish .	☐
■ Click Home to leave the configuration.	☐
■ Close the local service UI	☐



See also (→ Microsoft info, [https://technet.microsoft.com/en-us/library/cc786468\(v=ws.10\).aspx](https://technet.microsoft.com/en-us/library/cc786468(v=ws.10).aspx)) for password requirements and allowed characters.

4.1.6.2 Enable Autologin


- If on site with local access to the service software, and if password for the meduser account is known

ACTION	MRAWP/MRWP
■ Select Options ⇒ Service ⇒ Local Service .	☐
■ Enter the valid service key (service level 4 or 7) and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select Local Host Users	☐

ACTION	MRAWP/MRWP
Mark the Enable Autologin checkbox to enable the Autologin function. Fig. 289: SESO users 	☐
Click Save.	☐
Click Finish.	☐
Click Home to leave the configuration.	☐
Close the local service UI	☐

- If AutoLogin is disabled and meduser password is not known, use Service Software Utilities from Remote (Limited Permanent) to enable AutoLogin.
To do this, LAN or SRS connection must be possible (IP or host name is known, Service Key)

ACTION	MRAWP/MRWP
Remote login via SRS or LAN.	☐
Enter the valid service key (service level 3 or 7) and click OK.	☐
Select Utilities in the service platform.	☐

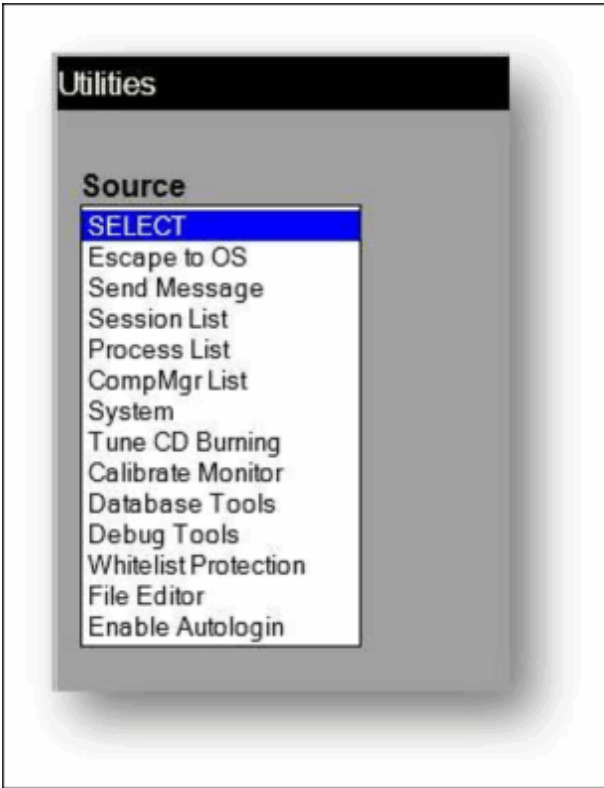ACTION	MRAWP/MRWP
On the Utilities page, under Source, select the entry Enable Autologin. <i>Fig. 290: Source entries</i> 	☐
Click Go. <i>Fig. 291: Utilities Autologin success</i> 	☐
Click Home.	☐
Close the service UI	☐


Now autologin is enabled and the password can be set again.

4.1.7 Advanced User handling

4.1.7.1 Enable Advanced User



The login for the advanced user can now be enabled or disabled. Also a password must be set if the advanced user is enabled in the Service Software UI.

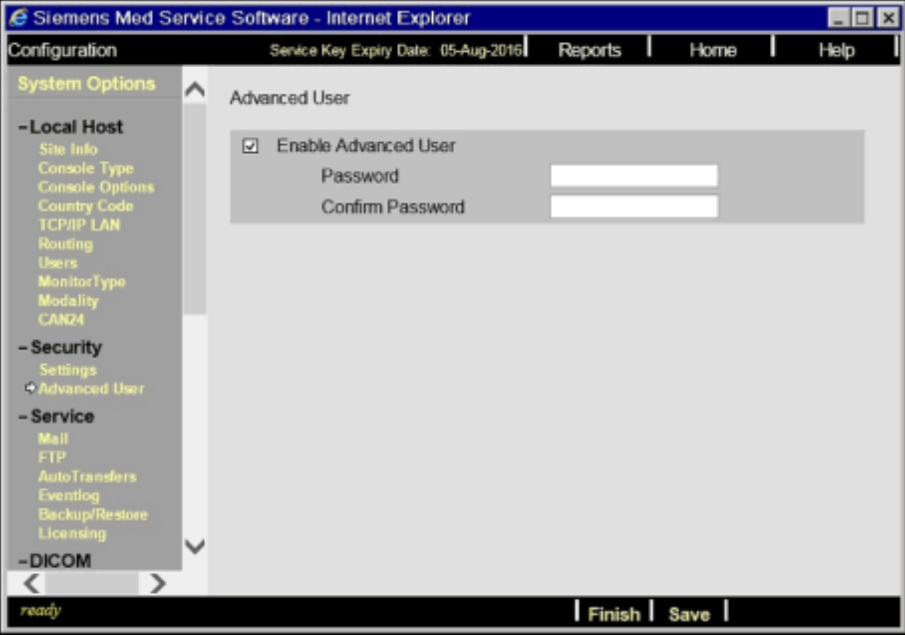ACTION	MRAWP/MRWP
■ Select Options ⇒ Service ⇒ Local Service .	☐
■ Enter the valid service key (service level 7) and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select Security Advanced User	☐
■ Mark the Enable Advanced User checkbox to enable the Advanced User login in the Windows Start menu. <i>Fig. 292: SESO Advanced User</i> 	☐
■ The password for the advanced user must now be set. The password must meet the following conditions: - Minimum length of 14 characters - It should have 2 uppercase letters and 2 lowercase letters - It should have 2 numbers and 2 special characters. The following 4 special characters are not allowed for the advanced user password: \ * ? ^	☐
■ Click Save .	☐
■ Click Finish .	☐

ACTION	MRAWP/MRWP
■ Click Home to leave the configuration.	☐
■ Close the local service UI	☐

4.1.7.2 Advanced User login

Procedure to display the complete Start menu:

- Run the Advanced User program from the Start menu (<CTRL> + <ESC>):
Start -> Programs -> Advanced User.

- Enter the user password:

You now have access to the complete Start menu.


4.1.7.3 Disable Advanced User



The login for the advanced user can now be enabled or disabled. No password must be set if the advanced user is disabled in the Service Software UI.

ACTION	MRAWP/MRWP
■ Select Options ⇒ Service ⇒ Local Service .	☐
■ Enter the valid service key (service level 7) and click OK .	☐
■ Select Configuration in the service platform.	☐
■ Select Security Advanced User	☐
■ Remove the mark from the checkbox Enable Advanced User to disable the Advanced User login in the Windows Start menu.	☐

Fig. 293: SESO Advanced User

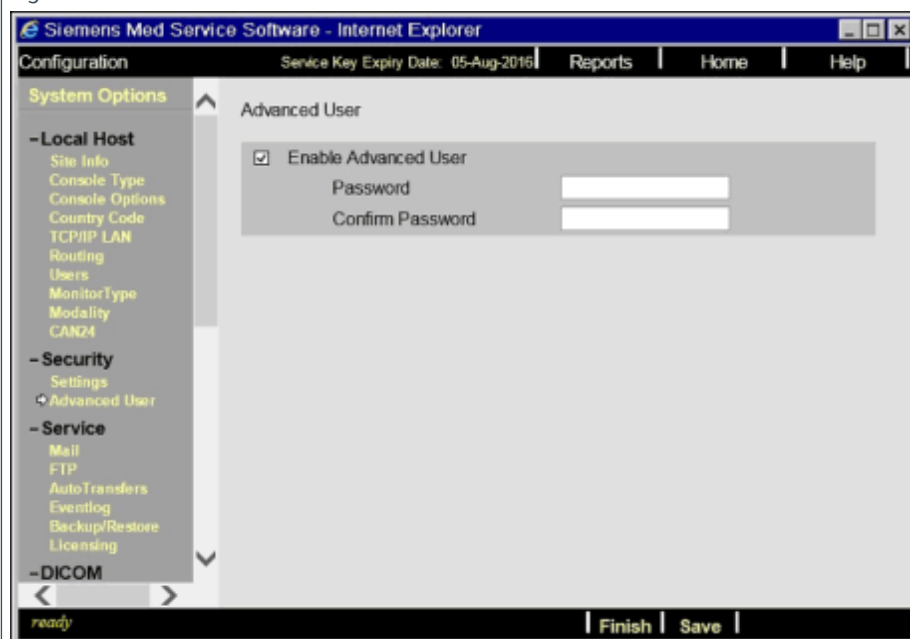


ACTION	MRAWP/MRWP
■ Click Save .	☐
■ Click Finish .	☐
■ Click Home to leave the configuration.	☐
■ Close the local service UI	☐

4.1.8 New feature for DICOM node configuration

With the new syngo version VH22B it is possible to configure the connection mode for each DICOM node (Network Nodes, Print Devices, HIS/RIS Nodes). If under connection mode the radio button secure (default) is selected, the DICOM communication is using the SSL/TSL protocol for the connection.

Fig. 294: DICOM Nodes configuration

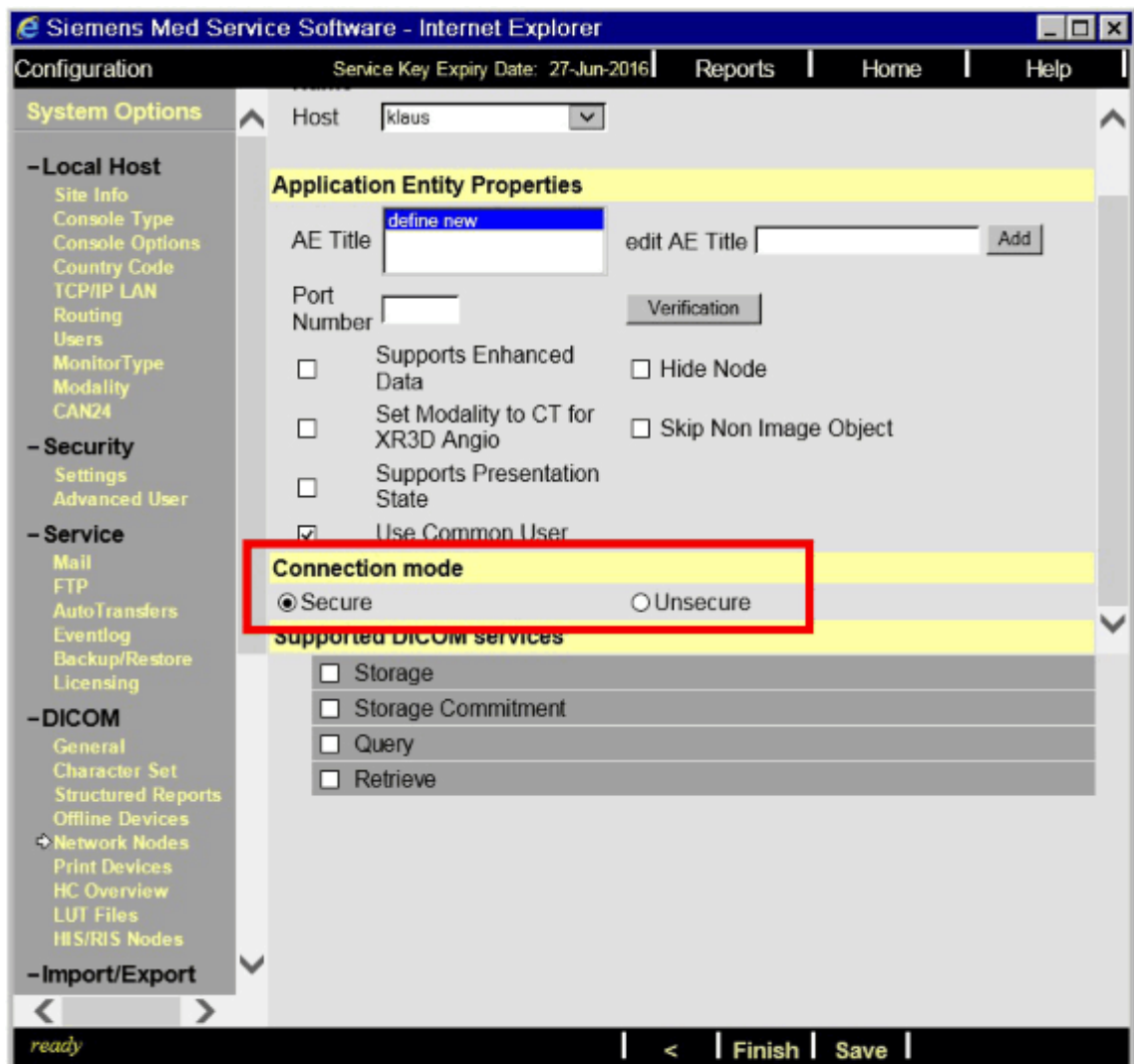

The connection fails if the DICOM node is not supporting the secure DICOM connection. In the event log the following message can be found (example).

151030|17:43:45|CSA_DSA|I16|TLS Handshake started with CONNECTOR. Proposing Cipher suite AES128-SHA:DES-CBC3-SHA


151030|17:43:45|CSA_DSA|E17|TLS Handshake rejected by CONNECTOR.

151030|17:43:45|CSA_DAL|E1|Error in Dicom Association Layer: An error occurred in Application 1001 within Association 0

In this case modify the connection mode from secure to unsecure for this DICOM node.

Please be aware that the default for the connection mode is secure if a new DICOM node is configured in the Service UI under configuration.

Fig. 295: Default

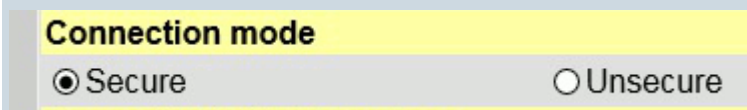

The connection mode configuration must be changed to unsecure. This is currently the general use case for the DICOM node configuration.

Use secure mode only if the DICOM node supports the secure DICOM communication with the SSL/TSL protocol.

4.1.9 System management installation workaround

4.1.9.1 Description:

If system management under Agent Controls is deinstalled, if the checkboxes

- SD -Asset - Management Agent
- Event Management Agent

are not marked (→ Fig. 296 Page 235), or if the installation is started again (both the checkboxes are marked), the Event Management Agent installation fails (→ Fig. 297 Page 236).

Fig. 296: MNP deinstalled


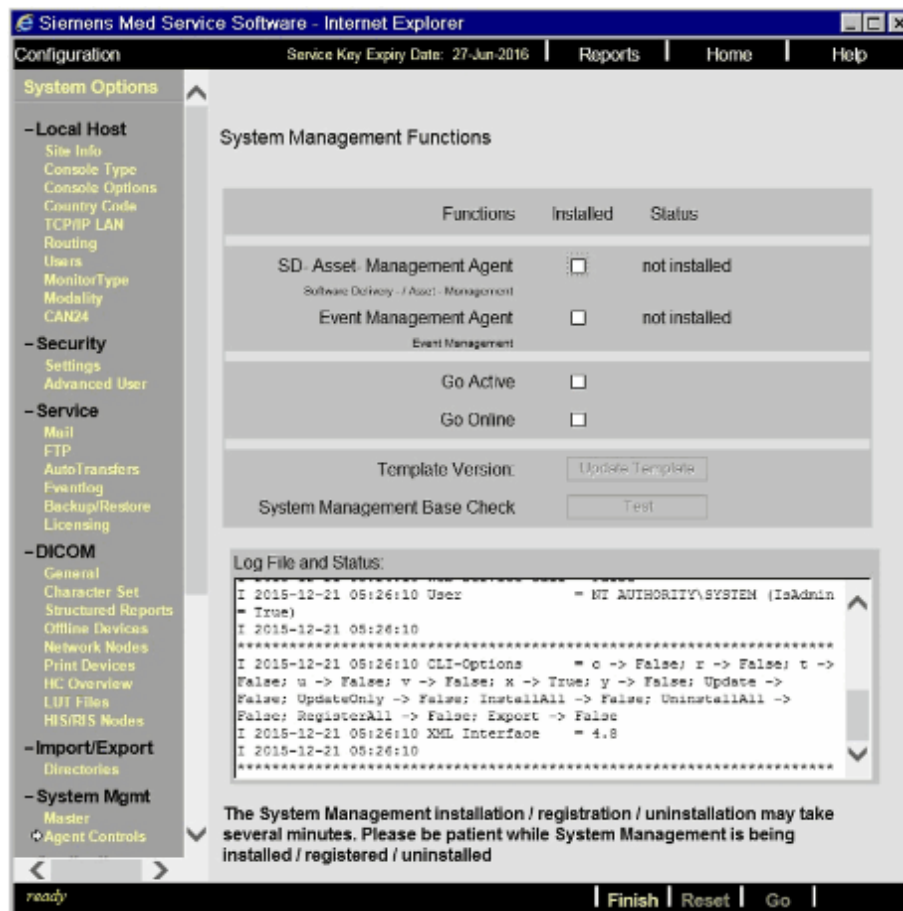


Fig. 297: MNP Event Agent installation failed



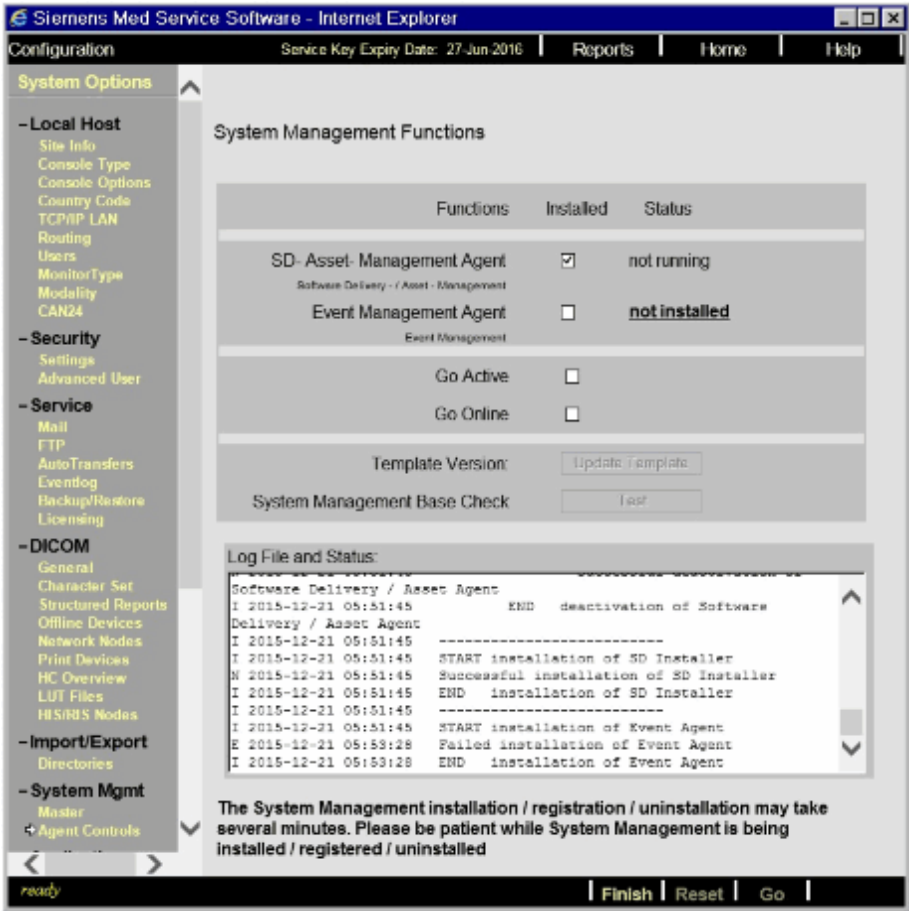


4.1.9.2 Solution:

The following procedure describes the workaround for this.

- Set Embedded Control in update mode

ACTION	MRAWP/MRWP
■ Close the local service UI.	☐
■ Enable Advanced User	☐
■ Start the MREmbeddedControlGUI from the Start menu Start ⇒ All Programs ⇒ MREmbeddedControl ⇒ MREmbeddedControlGUI	☐

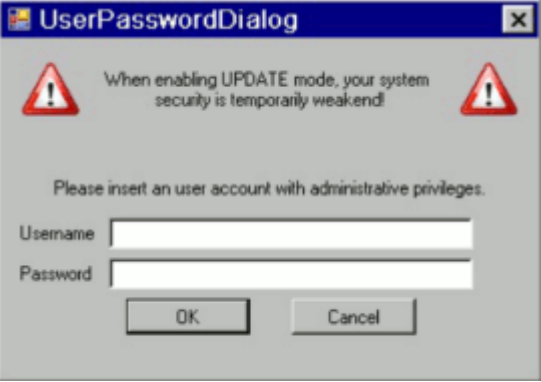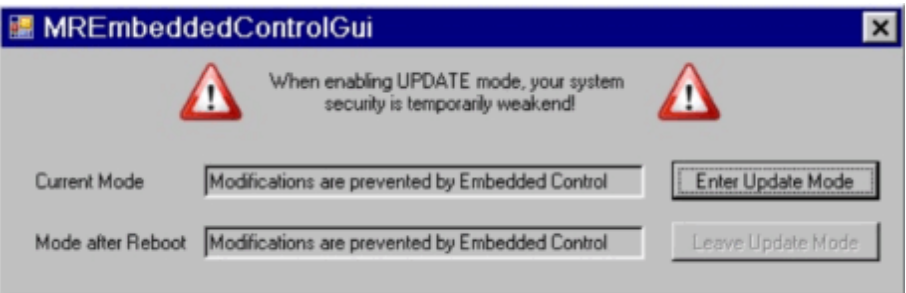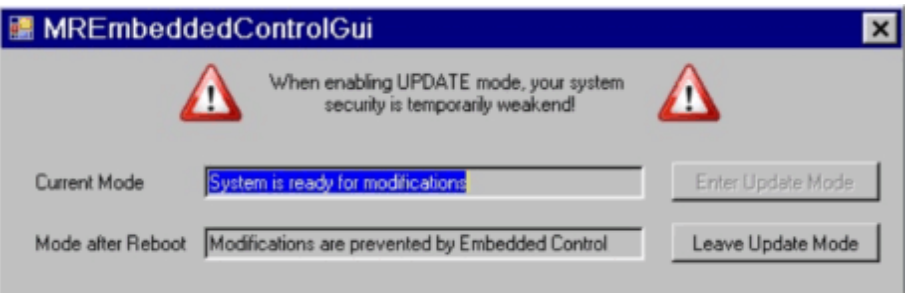
ACTION	MRAWP/MRWP
In the popup displayed, the user account with administrative privileges is requested. <i>Fig. 298: EC User Password request</i>  Type in as user "syadmin" and enter the password.	☐
Now the MREEmbeddedControlGUI is displayed <i>Fig. 299: MREEmbeddedControlGUI</i> 	☐
Press the button Enter Update Mode to switch embedded control into the update mode.	☐
Confirm the displayed popup with OK	☐
Now switch embedded control is in update mode. <i>Fig. 300: MREEmbeddedControlGUI update mode</i> 	☐
Leave the MREEmbeddedControlGUI open, don't close it.	☐

- System Management installation and activation

ACTION	MRAWP/MRWP
■ Select Options ⇒ Service ⇒ Local Service .	☐
■ Enter the valid service key.	☐
■ In the local service UI, select Configuration .	☐
■ Under System Mgmt. Agent Controls , set the following configuration: Mark the checkboxes "SD- Asset- Management Agent" and "Event Management Agent".	☐
■ Click Go .	☐

Fig. 301: MNP installation new

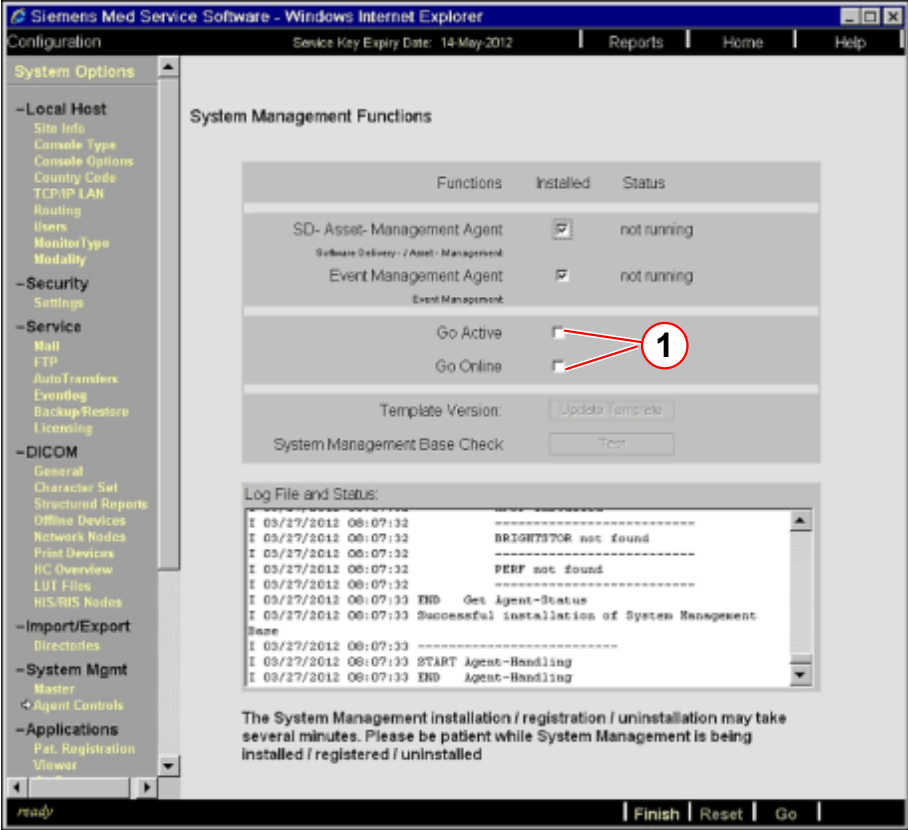


ACTION	MRAWP/MRWP
Now both Agents are installed and system management can be activated For the activation set the following configuration: Mark the checkboxes "Go Active" and "Go Online". (The other checkboxes have already been marked). Fig. 302: System Management activation  (1) Mark both checkboxes In the Log File and Status area, the current status of System Management activation/installation is displayed.	☐
Click Go.	☐
Note: The installation/registration will take approximately 10 to 15 minutes.	

ACTION

MRAWP/MRWP

- If activation is successful, the following is displayed:

Fig. 303: System Management activation OK

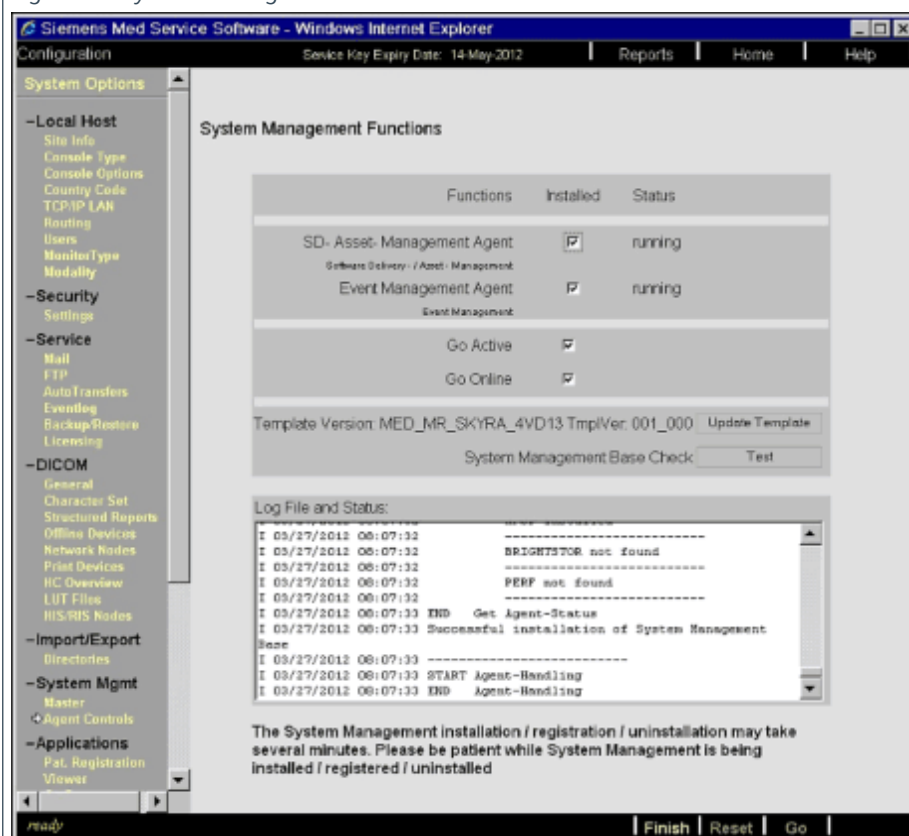


Select **Finish** to continue

- Click **Finish** and close the local service UI.


In the event of problems (installation/activation fails a second time), contact the SRS Help Desk (Tel. xx49-9191-188080-191). For help with problems within the US, call the US Uptime Service Center/Technical Support Group.

- Switch embedded control back into normal mode.

Fig. 304: MREmbeddedControlGui update mode

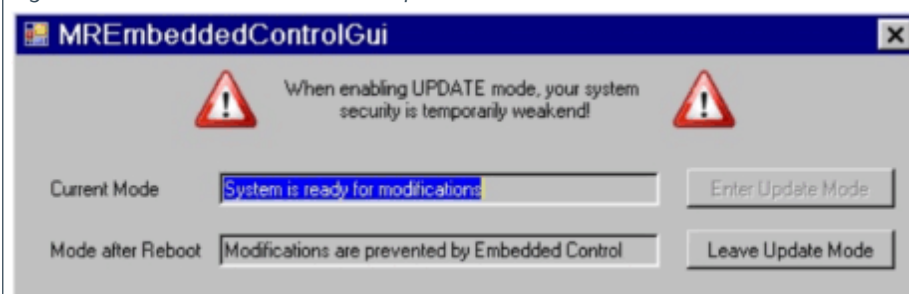

Press the button **Leave Update Mode** and confirm the displayed popup with **OK**

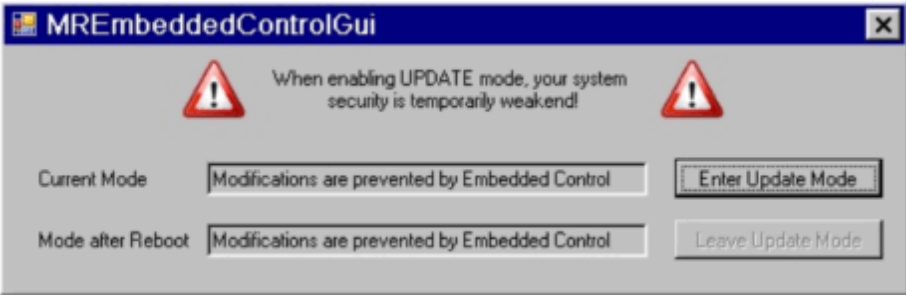ACTION	MRAWP/MRWP
■ Now EmbeddedControl is in normal mode. Fig. 305: MREmbeddedControlGUI 	☐
■ Close the MREmbeddedControlGUI by pressing the X	☐

4.1.10 Memory upgrade to 32 GByte



Use the ESD kit when performing the upgrade.

- Only valid for MRAWP HPz420 host with
- part number **10864002** (see Siemens part number sticker)
- or
- only valid for Brazil: part number **10928155** (see Siemens part number sticker)

Fig. 306: Example of partnumber sticker (Model No.)



4.1.10.1 Disassemble HPz420

- Remove all connected cables
- Use ESD outfit.

- Open the casing
 - Remove the two lock screws on the back of the casing.

Fig. 307: HP Z420 rear view




(1) Lock screw

(2) Lock screw


- Lay the device on its side on a stable, flat, and clean surface.
- Pull up on the handle (1), slide the cover 1/2 inch toward the rear of the system (2), then rotate the cover off the chassis (3).

Fig. 308: HP Z420 side panel

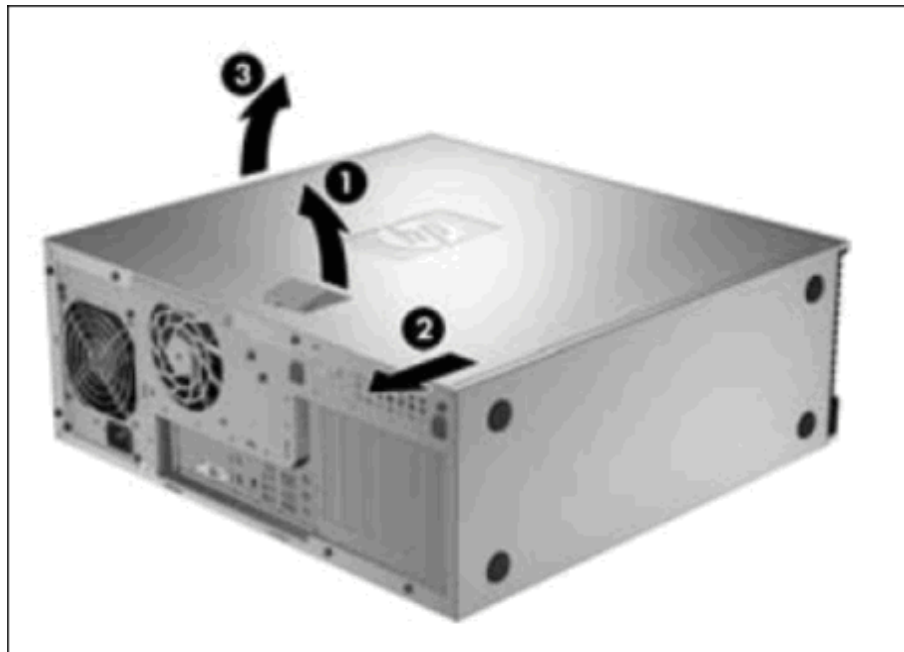

- Removing the fan from the second memory location


Fig. 309: HPz420 inside

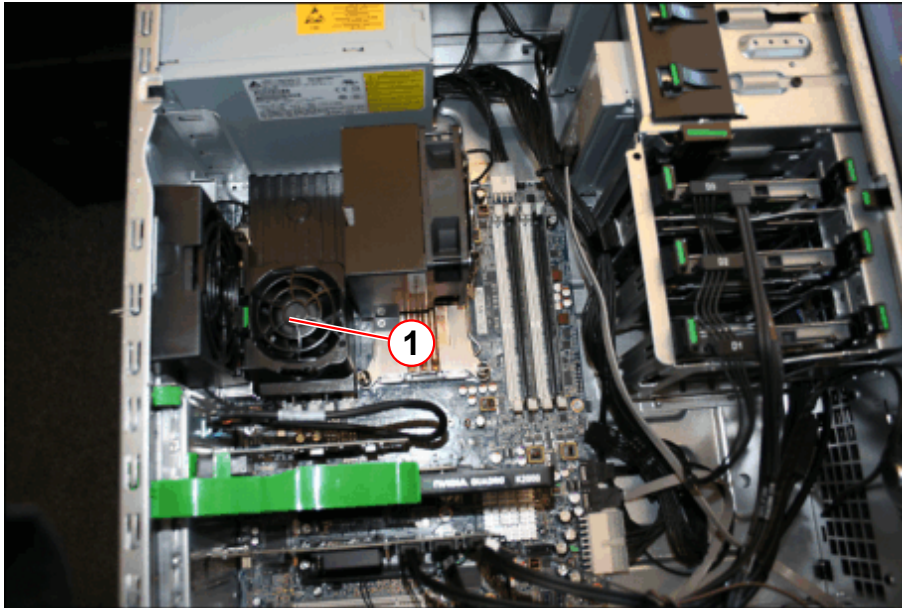

(1) Fan to remove


- Press the fan release latch as shown in the following figure, and guide fan out of the mounting location.

Fig. 310: Fan remove

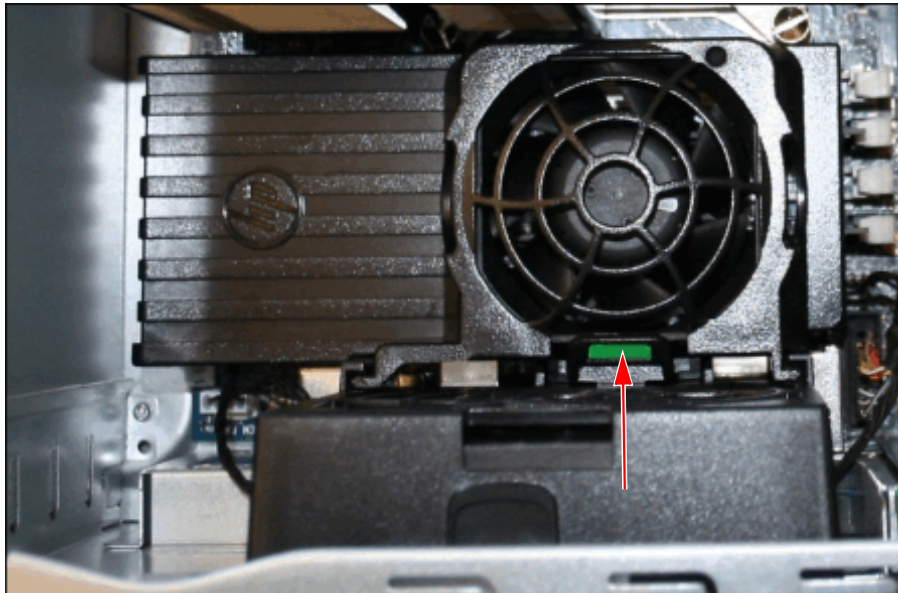

- Place the removed fan on the power supply (see following figure).

Fig. 311: HPz420 with removed Fan

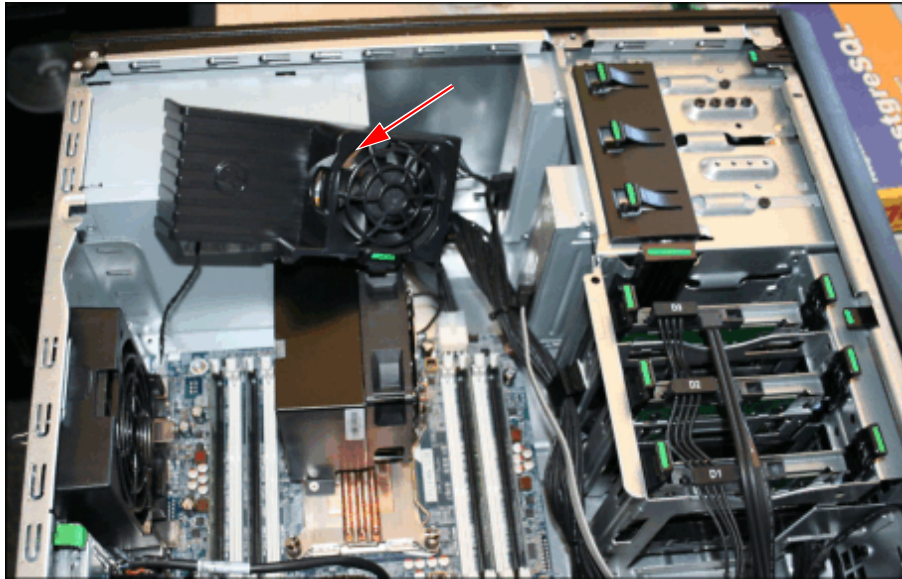

4.1.10.2 RAM installation for HPz420

- Install the additional 4 GB memory modules in the empty module locations (white sockets).




Only the following memory module can be used:

10928151 MEMORYMOD 4GB DDR3-1600 UNB (Siemens part number and description)

HP description: HP 4GB (1x4GB) DDR3-1600 ECC or DDR3-1866 UNB (both can be used)

Fig. 312: HPz420 memory module locations

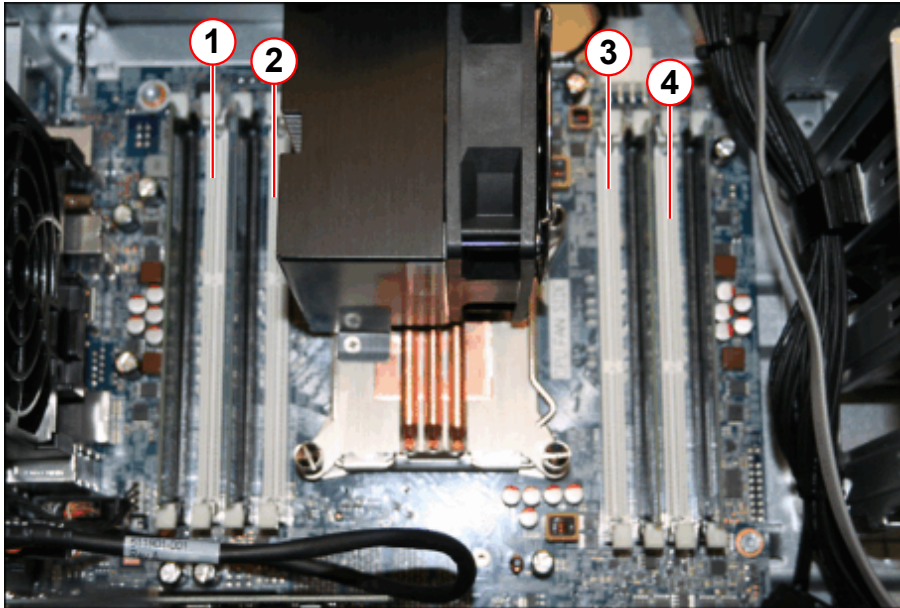

- (1) Place 4 GB memory module white empty DIMM socket
- (2) Place 4 GB memory module white empty DIMM socket
- (3) Place 4 GB memory module white empty DIMM socket
- (4) Place 4 GB memory module white empty DIMM socket

Insert the 4 GB memory modules in white empty DIMM socket.


When installing memory modules, make sure the modules align properly with the memory socket. There should be keys (small indents) on your memory modules that

fit the keys in the memory socket. DDR modules and sockets have only one key which is near the center of the module/socket.

Procedure to install a DIMM memory module

- Push gently outward on the DIMM socket levers as shown in the following illustration. (The example shows a HP Z600 Workstation.)

Fig. 313: Opening DIMM socket levers

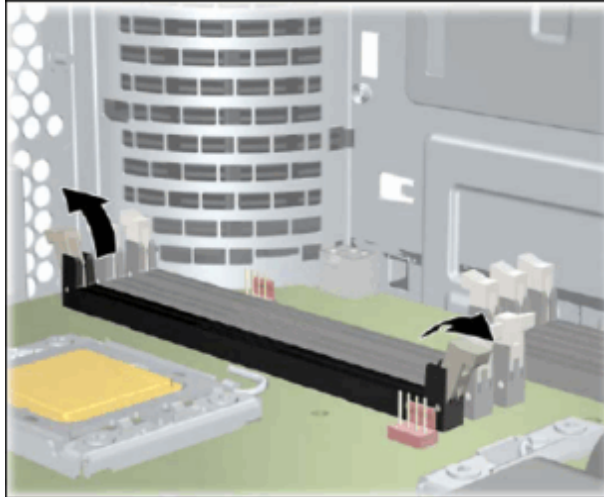


- Align the DIMM connector key with the DIMM socket key, and then seat the DIMM firmly in the socket (1) as shown in the following illustration.



DIMMs and their sockets are keyed for proper installation. To prevent socket or DIMM damage, align these guides properly when installing DIMMs.

Fig. 314: Seating the DIMM

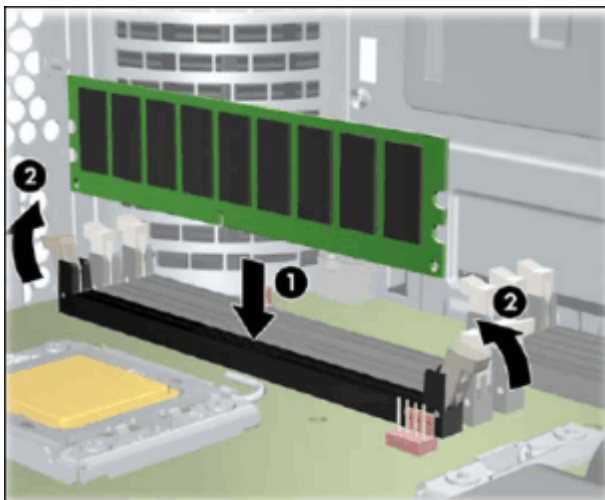

- Secure the socket levers (2).

4.1.10.3 Reassemble HPz420

- Installing the fan
 - Ensure that the latch on the fan snaps in.
- Closing the casing
 - Lay the device on its side on a stable, flat, and clean surface.
 - Hook the side cover into the cut-outs in the upper and lower guide rails. Pay attention to ensure that the lugs of the lock fit into the associated retainers on the casing.
 - Slide in the side cover until you hear it engage in the clip.
 - Secure the two lock screws of the side cover on the back of the casing.
 - Relabel the host HPz420 using the new part number (10928295, or for Brazil only, 10928294) sticker and new Rev. sticker.
 - Return the device to the place where you will need to use it.
 - Reconnect any disconnected cables (power cord, cables to external devices, etc.).


ACTION	MRAWP
■ Memory installed	☐
■ New part number label placed on host.	☐

4.1.10.4 Memory check

- Check recognized memory size
 - Ensure that the power cord connector is removed from the MARS.
 - Switch on the system.
 - Press F1 to skip the memory check from the BIOS.
 - Login as user administrator
 - On the desktop, select the icon "Computer"
 - Press the right mouse button and select *Properties* in the menu displayed.

- In the Properties window displayed, check the value for "Installed memory (RAM)" - it should be:
 - 32.0 GB

Fig. 315: System properties

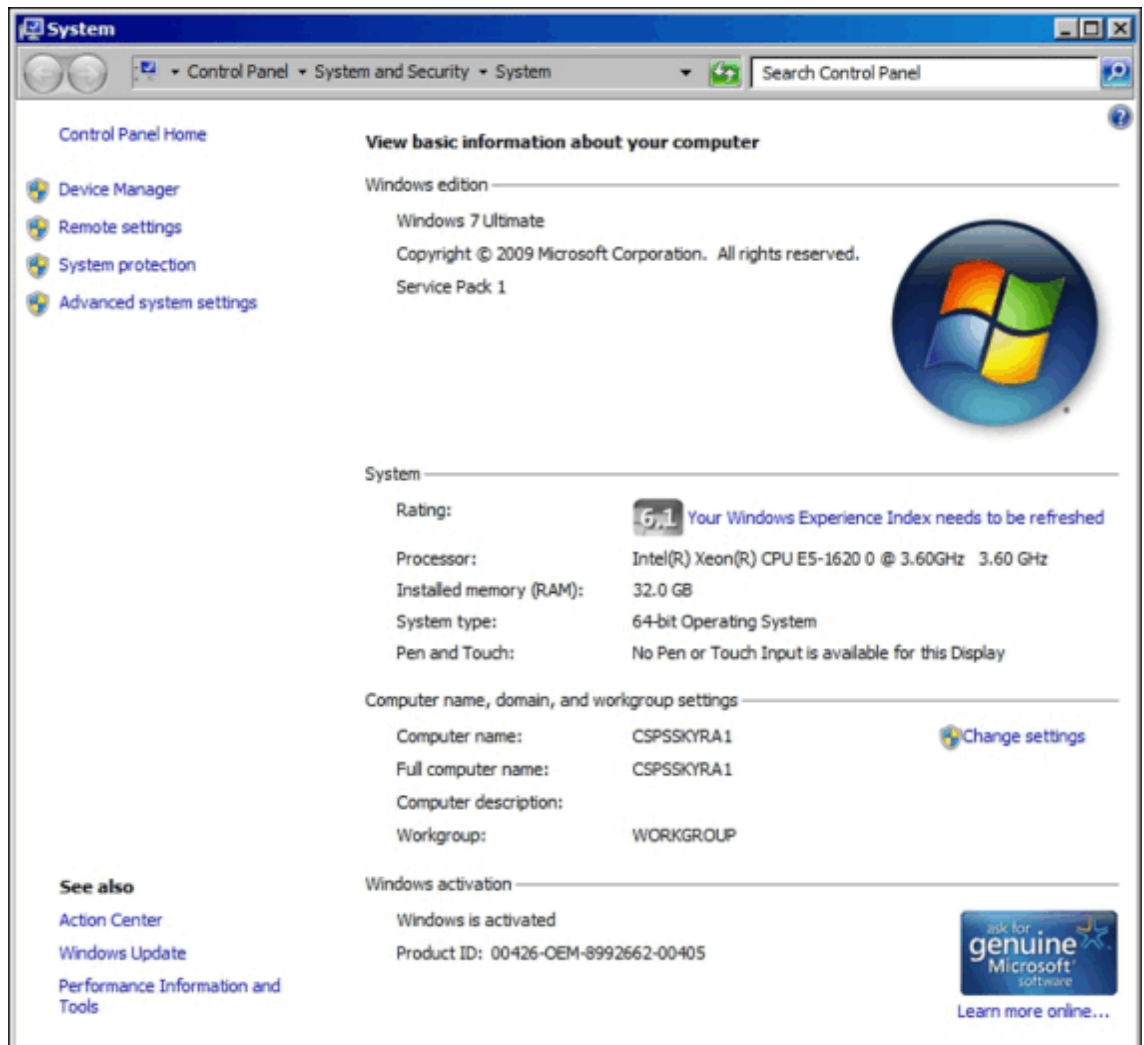

If this is not the case, the memory modules are not correctly installed.

- Close System Properties popup window.
- Continue with (→ MRAWP (MRC) BIOS Update / Page 121)

Version 02

- (→ Special Tools/Documents / Page 6) ESD kit added.
- (→ MRAWP (MRC) Memory upgrade to 32 GByte / Page 62) updated.
- (→ MRAWP (MRC) CSSI Installation / Page 71) workaround added.
- (→ Final Check - to be performed for MRC and MRSC (if applicable) / Page 110) Software label changed.

Version 03

- (→ Prerequisites / Page 5) added required UI's for disable OLM.
- Added update for software version VE11B

Version 04

- (→ Prerequisites / Page 5) added additional note concerning the OLM deactivation.
- (→ Content of the Update Kit / Page 7) Added the OLM material
- Added new chapter (→ Deactivation of OLM (if applicable) / Page 110)
- (→ Final Check - to be performed for MRC and MRSC (if applicable) / Page 110) Added new verification of the new version label

Version 05

- Added update for software version VE11C

There are no Hazard IDs in this document.

All documents may only be used for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

healthcare.siemens.com/services

Siemens Healthineers Headquarters
Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen
Germany
Telephone: +49 9131 84-0
siemens.com/healthineers

Print No.: MR00-000.814.56.05.02 | Replaces: MR00-000.814.56.04.02
Doc. Gen. Date: 02.16 | Language: English
© Siemens Healthcare GmbH, 2014

siemens.com/healthineers