

DQA CALIBRATION FOR 15X

PURPOSE:

In many HDxt systems running 15x the software does not work to run the automatic DQA calibration tool and it will not perform the scans to run the measurements. With this method we can run the scans manually and let the tool perform the measurements and calculate the changes.

1. Start a new exam and enter the following new patient data:
 - Patient id: geservice
 - Patient name: dqa
 - Weight: 111 lbs
2. Under the service protocols in other, choose adqa and click in *accept* to load all three series (for 3T you will need to run for WHOLE and ZOOM).




- Click on *save series* and *scan*.
- Click on *view edit* and *save series* in the following to series. Select *autoscan* to run them all automatically (it will save you from having to click *download* and *scan* after each series is done).
- When all three series are done, end the exam.


6. Navigate to the service tools and open the service browser if it is not open.
7. Select *DQA calibration* under the *Calibration* tab.
8. Click on *Click here to start this tool* to open up the tool.
9. When the tool opens click on *Start* at the bottom of the tool.
10. When the tool starts running, navigate back to the scan area and start a new patient.
11. Type *mr* in the patient id and it will reload all the values from before.
12. Under the service protocols select the *gradcal* protocol and select only the *dqa zfield cal*.



13. Click on *Save Series and Scan*.
14. Navigate back to the service tools and you should see the DQA tool running automatically.
15. Wait until it finishes all the calculations.

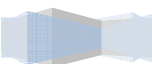

16. If all the calculations are within specification the tool will finish and it will show *completed* in green.

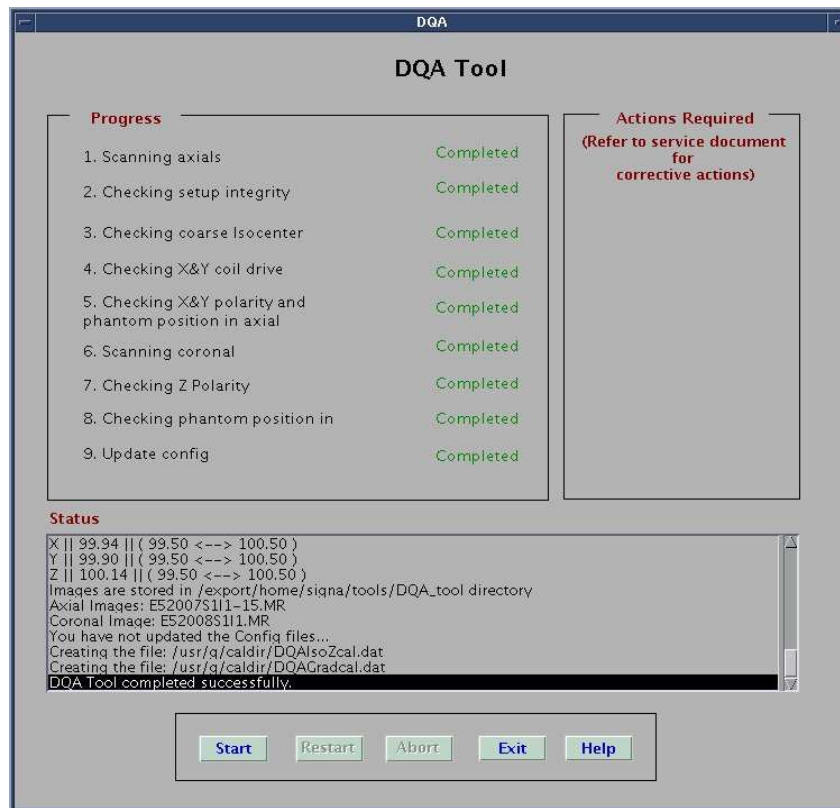

17. If any adjustment is needed, in most cases, the tool will fail. Write down the values and manually modify them in the *Config File Manager* (make sure to write down the original values).

18. Perform a TPS reset and rerun all scans and the tool again.

