

GE HealthCare

Technical Bulletin



Technical Bulletin Unique Identifier	DOC2953508
Title	Service Note of SIGNA MR355/MR360/MR380 new Mcat cables dummy option, new middle length gradient cable and available cable length clarification(China Region)
Subject	Service Note of SIGNA MR355/MR360/MR380 new Mcat cables dummy option, new middle length gradient cable and available cable length clarification(China Region)
Applicability (e.g., affected product, process)	China Region Only: 1.5T SIGNA MR355 1.5T SIGNA MR360 1.49T SIGNA MR380
Date	2-Feb-2024
Relevant Product Business (Select one)	☒ Imaging ☐ Patient Care Solutions ☐ Ultrasound ☐ Other: _____
Originating Modality	MR
Details	Service Note of SIGNA MR355/MR360/MR380 new Mcat cables dummy option, new middle length gradient cable and available cable length clarification(China Region)
Tracking Number	ECR2370861
Affected Service Publications	SIGNA MR355/ SIGNA MR360 Pre-installation manual (Class A: 5538859-1EN Rev3) SIGNA MR380 Pre-installation manual (Class A: 5538860-1EN Rev3)
Resolution	Service Note of SIGNA MR355/MR360/MR380 new Mcat cables dummy option, new middle length gradient cable and available cable length clarification(China Region)
Technical Bulletin Owner	Yuan Jing (SSO: 212476216)

(This SN is classified as Class A)

*** Printed copies are uncontrolled unless otherwise identified ***
 Before using this document, consult MyWorkshop for the latest revision.
 GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare

Technical Bulletin



1. Effectivity

China Region Only:

1.5T SIGNA MR355

1.5T SIGNA MR360

1.49T SIGNA MR380

2. Design overview:

Background:

According to China HPM service team request, MR355/360/380 shall have same strategy with 3T products), need set up dummy catalog into PCM&MCR which will help select the correct cable length after sitting measurements. Meanwhile help about reduce cable routing time and improve customer satisfaction. Creator/Explorer are working on the ECR to have same changes and available cable length clarification in couple weeks.

New Mcat dummy option and available cable length clarification


New/Legacy Mcat	Mcat Description	Available Cable Length mm(in.)	From	To	Length Identifier (shown in the figure below)
New	M50030SC Signa Creator Explorer Scan Room Kit	Dummy Option	ISC, bottom	MAG, bottom back	L3 (This length is for system cables)
Legacy	M50022SS Signa Creator Explorer Short Scan Room Kit	7800(307) Change from origin 4000 (157.44)			
Legacy	M50022SL Signa Creator Explorer LONG Scan Room Kit	11500(453) Change from origin 9200 (362.16)			

GE HealthCare

Technical Bulletin



New/Legacy Mcat	Mcat Description	Available Cable Length mm(in.)	From	to	Length Identifier (shown in the figure below)
New	M50030GC Gradient cables and Ground cables kit	Dummy Option	ISC, bottom	MAG, bottom back	GRADIENT CABLE(GND)
Legacy	M50022GS Gradient cables and Ground cables kit-short	4800(189)			
New	M50022GM Gradient cables and Ground cables kit-Mid	7800(307)			
Legacy	M50022GL Gradient cables and Ground cables kit-Long	11500(453)			

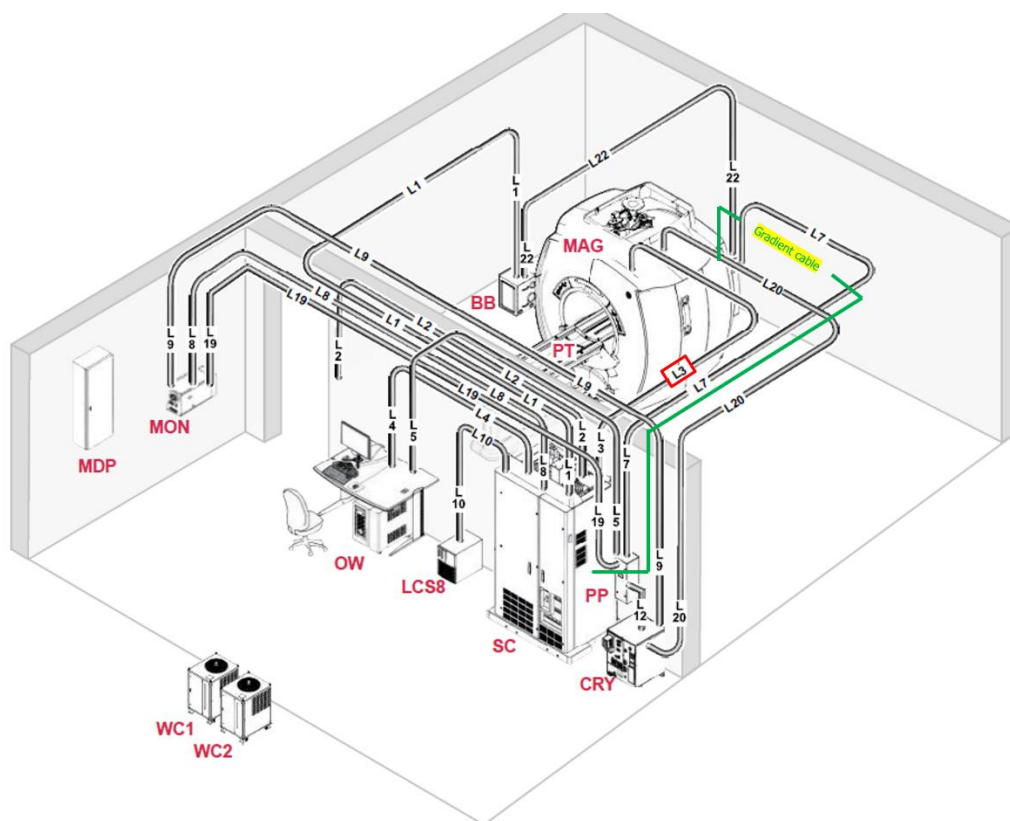

*** Printed copies are uncontrolled unless otherwise identified ***
 Before using this document, consult MyWorkshop for the latest revision.
 GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare

Technical Bulletin



MCR:

Must choose one set of Equipment Room and Scan Room cables.

Default option will be M50030SC and M50030GC. Selection must be updated during Cable Select Process.

Scan Room Kit	
CAT	Description
M50030SC	Signa Creator Explorer Scan Room Kit
M50022SL	Signa Creator Explorer LONG Scan Room Kit
M50022SS	Signa Creator Explorer Short Scan Room Kit

Gradient Cable Kit	
CAT	Description
M50030GC	Gradient cables and Ground cables kit
M50022GL	Gradient cables and Ground cables kit-Long
M50022GM	Gradient cables and Ground cables kit-Mid
M50022GS	Gradient cables and Ground cables kit-short

3. Appendix:

Part#	RUN#	Description	Cable Long/Short	From	To	Length (mm)	From: I/O to Component	TO: I/O to Component	Usable Length	Group	Note
6611002	M0503	Cable Run_M0503 MR2-F5-F6 TO MG6 J1	Long	ISC	Blower	15000	1500	500	13000	L1	
6611003	M0504	Cable Run_M0504 MR2-F7-F8 TO MG6 J2	Long	ISC	Blower	15000	1500	500	13000	L1	
5441863	M3522	Cable for LED power box and magnet enclosure LED light	Long	ISC	Rear Pedestal	15000	1500	500	13000	L3.2	
6611097	M3527	Cable Run_M3527 EMERGENCY OFF SW TO PP1-J18	Long	PP	E-off1	27000	700	500	25800	NOT SHOWN(L14)	
6611001	M4502	Cable Run_M4502 MR2-GND TO MG6-GND	Long	ISC	Blower	15000	1500	500	13000	L1	
5835766	P2500	Cable P2500 RRX J5 to ICE J32	Long	ISC	Magnet Side	17000	1000	1500	14500	L3.1	
6613100	M4505	Cable, magnet side kit to cabinet PP	Long	ISC	Rear Pedestal	19500	200	500	18800	L3.2	

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.

GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

Released

GE HealthCare

Technical Bulletin



6611003-1	M0504	Cable Run_M0504 MR2-F7-F8 TO MG6 J2	Short	ISC	Blower	8000	1500	500	6000	L1	
6613100-1	M4505	Cable, magnet side kit to cabinet PP	Short	ISC	Rear Pedestal	9500	200	500	8800	L3.2	
6611097-1	M3527	Cable Run_M3527 EMERGENCY OFF SW TO PP1-J18	Short	PP	E-off1	17000	700	500	15800	NOT SHOWN(L14)	Origin was 10000
6611001-1	M4502	Cable Run_M4502 MR2-GND TO MG6-GND	Short	ISC	Blower	8000	1500	500	6000	L1	
6611002-1	M0503	Cable Run_M0503 MR2-F5-F6 TO MG6 J1	Short	ISC	Blower	8000	1500	500	6000	L1	
5441863-1	M3522	Cable for LED power box and magnet enclosure LED light	Short	ISC	Rear Pedestal	9000	1500	500	7000	L3.2	Origin was 6000
5835766-1	P2500	Cable P2500 RRX J5 to ICE J32	Short	ISC	Magnet Side	12000	1000	1500	9500	L3.1	Origin was 9000
5801995-2	M3530	Starter, Cable Run_M3530 G-Filter-Y to Magnet-Y, 6m	Short	ISC	magnet	6000	1000	200	4800	GRADIENT CABLE(GND)	
6611050-11	M4500	Cable Run_M4500 MS1 TO RF-COM-GRND	Short	ISC	magnet	6000	1000	200	4800	GRADIENT CABLE(GND)	
5801996-2	M3531	Starter, Cable Run_M3531 G-Filter-Z to Magnet-Z, 6m	Short	ISC	magnet	6000	1000	200	4800	GRADIENT CABLE(GND)	
5801994-2	M3529	Starter, Cable Run_M3529 G-Filter-X to Magnet-X, 6m	Short	ISC	magnet	6000	1000	200	4800	GRADIENT CABLE(GND)	
6611050-10	M4500	Cable Run_M4500 MS1 TO RF-COM-GRND	Long	ISC	magnet	19500	1000	200	18300	GRADIENT CABLE(GND)	
5801996	M3531	Starter, Cable Run_M3531 G-Filter-Z to Magnet-Z	Long	ISC	magnet	12700	1000	200	11500	GRADIENT CABLE(GND)	
5801995	M3530	Starter, Cable Run_M3530 G-Filter-Y to Magnet-Y	Long	ISC	magnet	12700	1000	200	11500	GRADIENT CABLE(GND)	
5801994	M3529	Starter, Cable Run_M3529 G-Filter-X to Magnet-X	Long	ISC	magnet	12700	1000	200	11500	GRADIENT CABLE(GND)	

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.

GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

GE HealthCare

Technical Bulletin



5333375	M3514	Cable MG2-A42 J1-2 to MR2-J121	Common	ISC	Magnet side	20000	1500	2000	16500	L3.1
5333374	M3513	Cable MG2-A42 J1-1 to MR2-J120	Common	ISC	Magnet side	20000	1500	2000	16500	L3.1
5333353	M1515	Cable MG2-A41-J103 To MR2-J32	Common	ISC	Magnet side	18500	1500	2000	15000	L3.1
5333355	M1513	Cable MG2-A41-J106 To MR2-J29	Common	ISC	Magnet side	18500	1500	2000	15000	L3.1
5333354	M1514	Cable MG2-A41-J105 To MR2-J30	Common	ISC	Magnet side	18500	1500	2000	15000	L3.1
5333358	M3520	Cable MG2-A41-J11 To MR2-A11-J78	Common	ISC	Magnet side	19500	1500	2000	16000	L3.1
5333356	M3519	Cable MG2-A41-J5 To MR2-A11-J77	Common	ISC	Magnet side	19500	1500	2000	16000	L3.1
5393392	M3540	Cable from MG2-A41-J10 to MR2-A11-J79	Common	ISC	Magnet side	19500	1500	2000	16000	L3.1
5333387	M3528	MG3-A2-J6 TO PP1-J12	Common	ISC	Rear Pedestal	15000	1500	1000	12500	L3.2
5333386	M3534	MR2-A11-J74 TO MG3-A3-J6	Common	ISC	Rear Pedestal	15000	1500	1000	12500	L3.2
5333385	M3533	MR2-A11-J76 TO MG2-A12-A2-J6	Common	ISC	Magnet side	16000	1500	2000	12500	L3.1
5333384	M3532	MR2-A11-J75 TO MG2-A12-A2-J5	Common	ISC	Magnet side	16000	1500	2000	12500	L3.1
5333388	M3535	MR2-J11 TO MG3-A7	Common	ISC	Rear Pedestal	15000	1500	1000	12500	L3.2
46-328000G974	828	Run 828 MS1-A3-A1-J403 TO PP1-J10	Common	PP	magnet	18440	700	5000	12740	L7.1
6611096	E3500	Cable Run_E3500 FACILITY DISCONNECT TO PP1-J18	Common	PP	E-off2	19812	700	500	18612	NOT SHOWN(L13)
5390678-4	1501	RUN 1501 PP1-J17 TO OM1-J16	Common	PP	Oxygen Monitor	30632	700	500	29432	NOT SHOWN(L6)
5390678-5	1502	RUN 1502 PP1-J17 TO OM3-J4	Common	PP	Oxygen Monitor2	27400	700	500	26200	NOT SHOWN(L11)
5333387-2	M3528	MG3-A2-J6 TO PP1-J12 no-magnet connector	Common	PP	Rear Pedestal	15000	700	1000	13300	L7.2
5803282	E3503	Cable Run_E085 PP1-J12 TO OW1-A21-J7	Common	PP	OW	30000	700	500	28800	L5
5534139	825	Run 831 FJ1 TO MSM1-A1-J7	Common	PP	MON	24384	700	500	23184	L19

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.

GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.

Released

GE HealthCare

Technical Bulletin



5508123	824	Run 824TR PP1-J10 TO FJ1	Common	PP	CRY	24384	700	500	23184	L12
2155316-2	623	623 CBL COMPRESSOR TO PEN PNL, OBO 49.21 FT OR 15.00 M LONG	Common	PP	CRY	15000	700	500	13800	L12
2218292	624	624 CONTROL CABLE, COLDHEAD	Common	PP	magnet	15000	700	4500	9800	L7.1
6611010	E3506	Cable RUN_E3506 MR2-J15 TO RF-Door- switch	Common	ISC	Door- sensor switch	30480	500	500	29480	L2
5479529	E3501	Run E3501 ISC-IO-J26 to Run E3047, 30000mm	Common	ISC	OW	30000	500	500	29000	L4
5502144	E3042	RUN E3042 ISC-Cabinet IF-J112 to E3049, 30,000mm	Common	ISC	OW	30000	500	500	29000	L4
5508128	827	Run 827 FJ3 TO FJ4 and MSSA5	Common	MON	CRY	23220	500	500	22220	L9
5508124	826	Run 826 FJ3 TO MSM1- A1-J9	Common	MON	CRY	18440	500	500	17440	L9
5508125	823	Run 823 MSM1-J10 TO MR2-A11-J24	Common	ISC	MON	18288	500	500	17288	L8
5145249	2009	Cable Run_2009 MR2-J47 TO MG2-A29-J1	Common	ISC	magnet	16916	1500	2000	13416	GRADIENT CABLE(GND)
5780443	E3037	USB Extender Cable, 30000mm, Run E3037	Common	ISC	OW	30000	500	500	29000	L4
2128699	606	606 MAGNET END TO ERU MRU	Common	Magnet	MRU	30480	500	500	29480	NOT SHOWN(L18)
6611505	E3505	Cable Run_E3505 WC1 and WC2 TO MR2-J301	Common	ISC	LCS/MCS	10000	500	500	9000	L10
6611041	E0501	Cable Run_E0501 WC2 TO MR2-T2	Common	ISC	LCS/MCS	10000	500	500	9000	L10

Revision	Version	Date	Author	Description
01	01	2-Feb-2024	Yuan Jing (SSO: 212476216)	Initial Release

*** Printed copies are uncontrolled unless otherwise identified ***

Before using this document, consult MyWorkshop for the latest revision.

GE HealthCare Confidential

General service documentation.

© 2023 GE HealthCare – GE is a trademark of General Electric Company used under trademark license.

Reproduction and/or distribution is prohibited.