


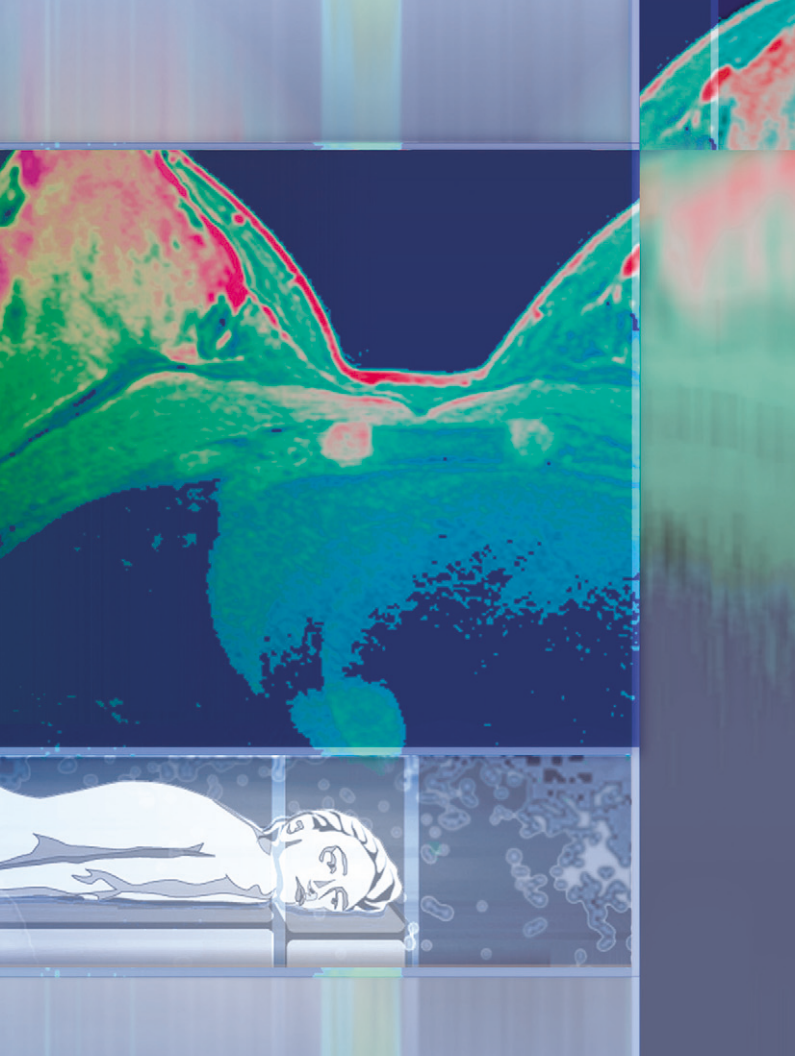

MR Mammography

Application Brochure

MAGNETOM Avanto
MAGNETOM Espree
MAGNETOM Symphony a Tim System
MAGNETOM Trio a Tim System
MAGNETOM Verio

www.siemens.com/healthcare

SIEMENS

© Siemens AG 2008–2009
All rights reserved

Siemens AG
Wittelsbacherplatz 2
D-80333 Muenchen
Germany

MR Mammography

Application Brochure

MAGNETOM Avanto

MAGNETOM Espree

MAGNETOM Symphony a Tim System

MAGNETOM Trio a Tim System

MAGNETOM Verio

This brochure describes MR mammography used together with syngo MR. It has been developed for medical personnel working in the field of MR tomography.

To optimize the user-friendliness of this brochure, its contents have been divided into four subject matters.

The first part of the brochure focuses on the basics and provides know-how about the subject matter.

The second part describes practical applications and their use based on sample examinations. The third and fourth part provides the basics as well as use with prosthetic silicone implants.

MAGNETOM Espree – Pink/*syngo* BreVis:

For the MAGNETOM Espree – Pink and the *syngo* BreVis application, several enhanced post-processing tools are available that are described in the BreVis manual.

Furthermore, the *syngo* BreVis software provides its own motion correction, subtraction, and time curve evaluation tools.

Contents

Basics	Overview of MR Mammography	2
MR Mammo- graphy	Sensitivity and specificity for breast lesions	4
	Qualitative evaluations for breast lesions	6
	Motion correction for dynamic contrast-supported measurements	8
	Quantitative evaluations for breast lesions	10
	Parameter images (Washin/Washout)	12
	Parameter images (TTP, PEI)	14
	Parameter images (MIPT)	16
Application	Dynamic examination overview	18
MR Mammo- graphy	Prepare contrast agent injection	20
	Diffusion-weighted measurements with REVEAL (optional)	22
	Set dynamic measurement (inline)	24
	Select parameter images (inline)	26

	Perform dynamic measurement (inline)	28
	Performing motion correction (if required)	30
	Reconstructing combination images (offline)	32
	Evaluate contrast agent dynamics (offline)	34
	Manual post-processing (offline)	36
	Manually scale subtraction images (offline)	38
Basics	Prosthetic silicone implants	40
Silicone implants	Spectral differences	42
Application	Displaying pure silicone images	44
Silicone implants	Perform manual adjustments (I)	46
	Perform manual adjustments (II)	48
	Index	50

Overview of MR Mammography

MR Mammography is the most sensitive method for detecting and staging invasive breast carcinomas. A special application is used in areas around prosthetic implants.

Suited for

- Verification of invasive carcinoma
- Prosthetic implants
(complications, exclusion of tumors)
- Evaluation of benignity and malignancy
in cases of structural thickening
- Preoperative determination of
multicentricity or multifocality
for breast-preserving therapy
- Distinction between mastopathies
and malignant tumor formations
- Differentiation of scar tissue
- Follow-up control after tumor surgery

Not suited for

- Evaluation of microcalcification
- Inflammatory changes
- Mammary secretions
- Known mastopathies

*Typical
breast
diagnosis
sequence*

**Medical history and
clinical examination**



Mammography



Sonography



MR Mammography

Application

MR Mammography is a valuable supplemental method for detecting breast lesions. MR Mammography is the method of choice when diagnosing prosthetic complications.

Sensitivity and specificity for breast lesions

To date, there is no general agreement regarding the role of MR imaging in the detection and characterization of breast lesions.

Sensitivity

Probability that the disease exists if the test results are positive.

Author	Sensitivity
Allgayer 1993	90%
Harms 1993	93%
Heywang 1989	100%
Allgayer 1995	96%
Kaiser 1993	99%
Boetes 1994	95%
Fischer 1993	95%
Kuhl 1999	91%
Drew 1999	100%

Extremely high sensitivity is achieved when the carcinoma measures 3 mm or more.

Specificity Probability that the test results will be negative, if the disease does *not* exist.


Author	Specif.	Initial	Curve	Morph
Allgayer 1993	29 %	–	–	–
Harms 1993	37 %	–	–	–
Heywang 1989	42 %	+	–	–
Allgayer 1995	51 %	–	–	–
Kaiser 1993	97 %	+	+	–
Boetes 1994	86 %	+	–	+
Fischer 1993	90 %	+	+	+
Kuhl 1999	37 %	+	–	–
Drew 1993	83 %	+	+	–

The specificity of MR Mammography is only moderate compared to the sensitivity. For this reason, limit contrast-enhanced MR imaging to specific indications and use it as an auxiliary method only.

Qualitative evaluations for breast lesions

For the detection of breast lesions, MR images as such are *not* informative. An evaluation of the dynamic image data is required to diagnose breast lesions.


*Subtraction
image*

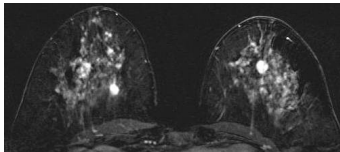

Subtraction

In subtraction images, areas enriched with contrast agent are displayed bright. In this case, a native data set (without contrast agent) is subtracted from a data set after administration of contrast agent.

MIP image

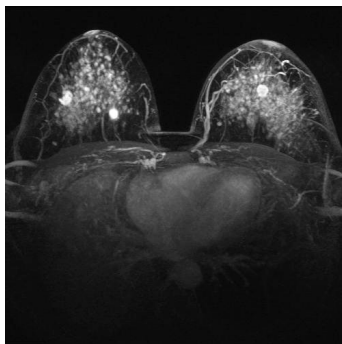

MIP

You can also calculate MIP images (MIP = Maximum Intensity Projection) to display vascular courses. The pixels with the maximum intensity are displayed in the MIP image.


Motion correction for dynamic contrast-supported measurements

To diagnose tumors in contrast agent-supported MR images, it has to be ensured that the patient does not move during the entire measurement.

Measurements with contrast agent usually take between 6–8 minutes. This means that motion artifacts may occur. These movements especially complicate diagnoses of very small (multifocal) tumors.

Motion correction

These artifacts can be corrected retroactively with motion correction. This means that even heavily blurred MR measurements can be used for the diagnosis.

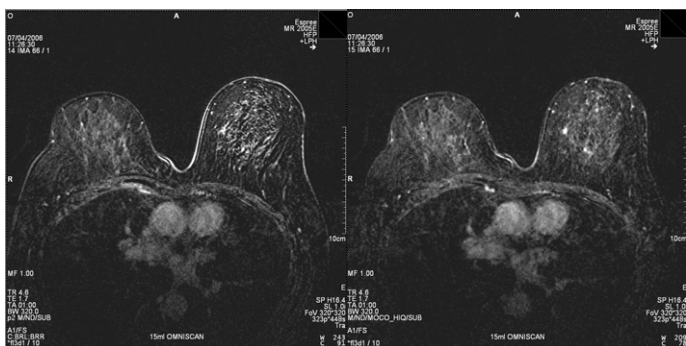

Subtraction analysis of a contrast-supported dynamic measurement:

without motion correction (left)

after motion correction (right)

Quantitative evaluations for breast lesions

Certain evaluation of the benignity or malignancy of tumors requires an additional examination of the temporal course of contrast agent enhancement.

Parameter images


Images that display the temporal course of contrast agent dynamics in color. Standard protocols from Siemens are optimized to spatial and temporal resolution and serve as the basis for parameter image reconstruction. In parameter images, pixels of different tissue types show a different signal behavior.

Parameters to characterize dynamics

- TTP: Time to maximum enrichment (*time to peak*)
- PEI: Area under the enrichment curve (*positive enhancement integral*)
- Washin: Signal increase in the area of enrichment
- Washout: Signal behavior after enrichment
- MIPt: Greatest signal increase during the dynamic measurement

**Combination
images**


Images that simultaneously display multiple parameters, typically Washin combined with Washout.

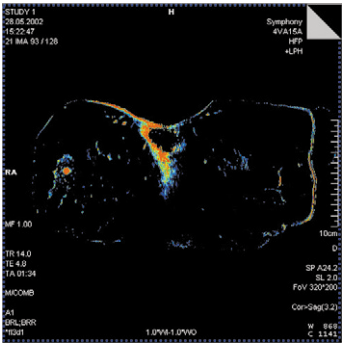

Washin/Washout combination image

Parameter images (Washin/Washout)


The Washin/Washout parameter images display signal changes at the starting and end areas of the dynamic curves. After contrast agent injection, breast carcinomas display a typical increase in signal intensity with earlier plateau formation and a wash-out effect.

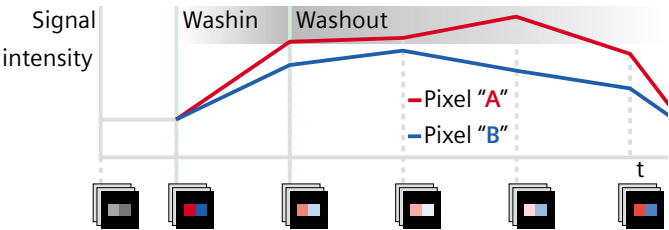

Pixel-by-pixel, the percentage change in signal value per minute is calculated based on the initial series, and is converted into gray and color values.

Highest value During the washin calculation, the parameter "Highest value" ensures that the signal change is measured up to the maximum. This prevents "false" results in case there is a series between the start and end series with a higher signal intensity than the end series.

Down sloping washout curve. Strongly negative values in the card are suspect. For this reason, a color palette with an inverted course of colors is used for the evaluation.

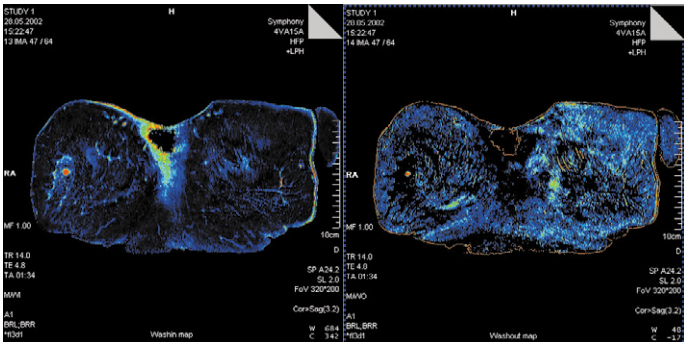


Washin image (left): The rise in signal in the starting area of the dynamic curve is displayed.


Washout image (right): The signal behavior at the end area of the dynamic curve is displayed.

Parameter images (TTP, PEI)


The time-to-peak image displays (pixel-by-pixel) the time in seconds until signal maximum. In a typical dynamic MR breast examination, multiple identical measurements (e.g., 6–8) are acquired, yielding just as many time points (6–8) for the dynamic curve evaluation and color maps.

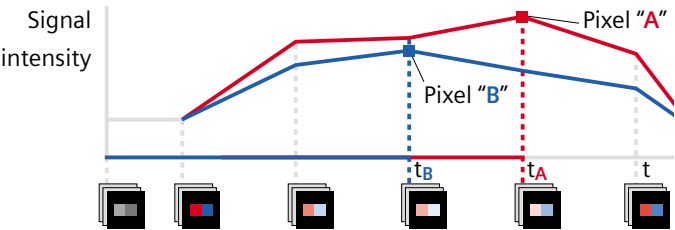

TTP image

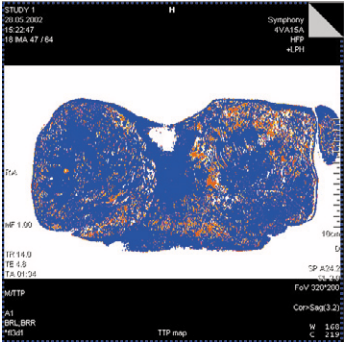


The positive enhancement integral image displays (pixel-by-pixel) the area under the time intensity curve.

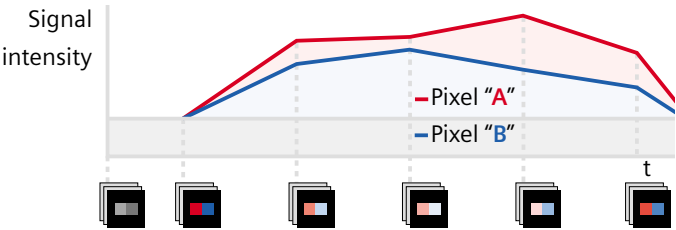


*Positive enhancement integral for two different pixels.
The value in the results image is greater for Pixel "A" than for Pixel "B".*

PEI image

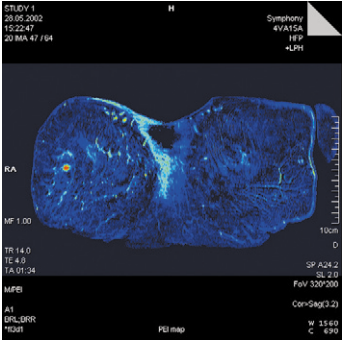


Parameter images (MIPT)

The MIPT calculates a Maximum Intensity Projection over time for each pixel. In the MIPT-image, you can see the extent of the increase in signal.

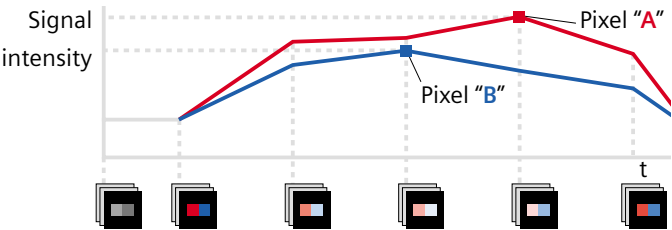


The brightest pixel in the dynamic time range is used for the results image.

MIPT image

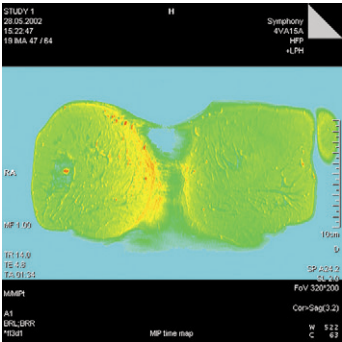

Dynamic examination overview

Precise detection of breast lesions is the goal.
For this reason, we measure T1-weighted precontrast and postcontrast images and then perform an evaluation (subtraction, parameter images, etc.).

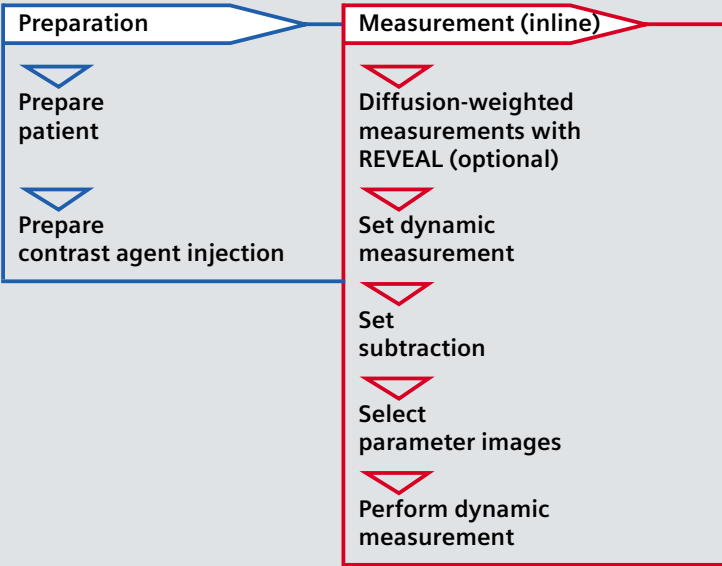

PROTOCOL

Localizer

3D FLASH dynamic

VIEWS 1+5M

(1 precontrast, 5 postcontrast measurements, incl. inline reconstruction of parameter images)

Post-processing



Perform motion correction (if required)



Reconstruct combination images



Mean Curve



Manual post-processing



Manually scale subtraction images

You can also use transverse slices with the GRAPPA iPAT option (PAT factor 2).

This ensures the same temporal resolution as with coronal slices. Additionally, transverse slices enable good evaluation of the ductal course.

Prepare contrast agent injection

Position the patient prone in the breast coil.

Recommendations. For MAGNETOM Espree and Verio: Position the breast coil on the spine matrix so that the breast is closer to the isocenter of the magnet.

Route the tube for the infusion *prior to* moving the patient table into the magnet.

Recommendations. The tube should be long enough so that it can be accessed from outside when the patient is in the magnet bore.


Insert a cannula into the forearm vein. Connect it to the injection tube with the attached 3 -way stop valve. Fill the contrast agent and saline solution in two separate syringes and connect the syringes to the stop valve outside the magnet bore.

Preparation

Prepare patient

Prepare contrast agent injection

Recommendations. Do *not* use I.V. bottles.

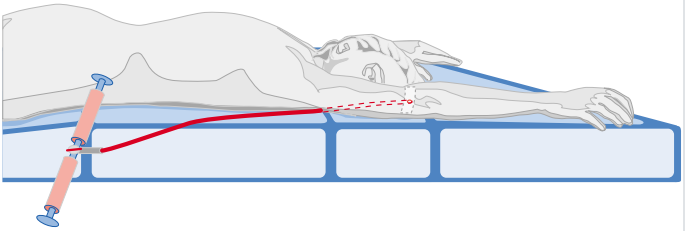

Instruct patient. Prior to moving the patient into the magnet, explain that she should remain completely still during the entire examination.

Diffusion-weighted measurements with REVEAL (optional)

For additional evaluation of lesions or follow-up, you are able to measure diffusion-weighted images with REVEAL. For this purpose, the lesion made visible in the dynamic examination has to be located in the diffusion-weighted images.

Recommendations. Since the contrast agent affects the T1 relaxation time and thus the effectivity of the fat suppression, an EPI-STIR measurement should be taken *before* the dynamic measurement.

REVEAL protocols. In the protocol tree under breast/library/other; with inline-ADC computation (prerequisite: Inline Diffusion option).

Determining the ADC value

The ADC map (ADC = Apparent Diffusion Coefficient) allows for an in-vivo pseudo quantification of the diffusion effect. Due to the limited mobility of water molecules in the vicinity of a tumor, malignant lesions show lower ADC values than benign lesions.

DWI with REVEAL




Transfer an ROI within the lesion from the diffusion-weighted images to the ADC map.

For this purpose, go to the Viewing task card and draw an ROI into the diffusion-weighted image.

Copy the ROI to the ADC map.
Determine the ADC value (mean value of the grey-scale values).

TIP

Drawing ROIs: Central necroses should not be enclosed by an ROI. It has a different cellularity and therefore a different ADC value compared to the tumor.

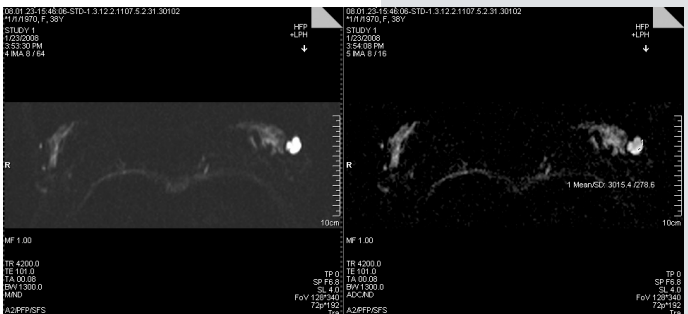


Set dynamic measurement (inline)

Localizer is measured/
dynamic
protocol is
open

Set any additional parameters desired in
the dynamic protocol.

Please note: Use multiple measurements
to prevent repeated adjustments between
measurements.

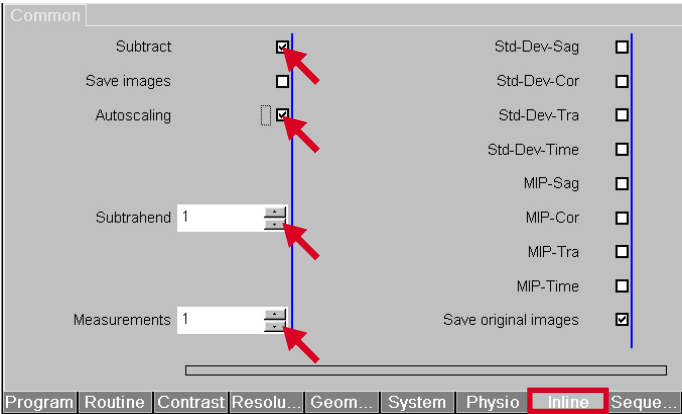

Inline card



Set dynamic
measurement
Set subtraction



Parameters for subtraction

The parameters are optimized in the Siemens standard protocol and usually may not be changed.

The number of the last measurement without contrast agent is set as the “Subtrahend”


“Autoscaling” ensures that any negative values that occur during the subtraction are *not* set to zero.

TIP

In addition, you can create MIP images in reconstruction directions sag, cor, and tra.

Select parameter images (inline)

If you would like to activate inline computation for dynamic protocols, please use the Inline Breast card.

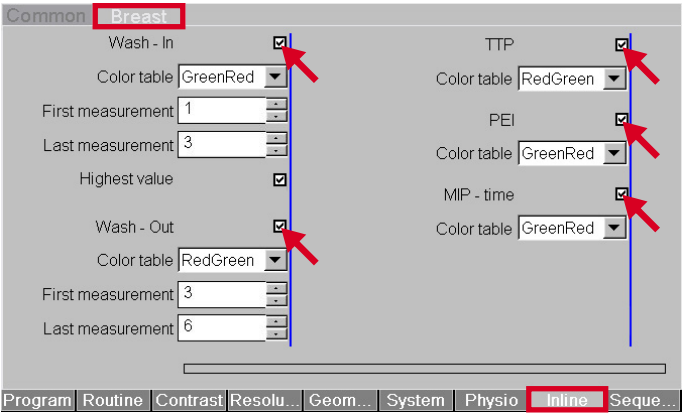

Inline Breast card

Select the parameter images desired using the Inline Breast parameter card.



Select parameter images



Reconstruct Washin image

Set the start value to the last precontrast series and the end value to the series in which the enrichment phase ended.

Reconstruct Washout image

Set the end value from the Washin parameter as the start value, and the last series as the end value.

TIP

You can also have the parameter images computed retroactively.

TIP

Windowing Washout images: Orientation to the contrast agent washout at the heart can be helpful.

TIP

Windowing can be made easier by using the **Perfusion** color table.

Perform dynamic measurement (inline)

Parameters for the dynamic measurement have been set

Check the infusion site and patency of the infusion tube immediately prior to starting the dynamic series.

Start the dynamic precontrast series.


Inject the contrast agent bolus immediately after the first measurement and then irrigate with a sufficient quantity of saline solution.

The subsequent postcontrast measurements start automatically.


You obtain dynamic precontrast and postcontrast series as well as subtraction and parameter images.

Measurement (inline)

Perform dynamic
measurement

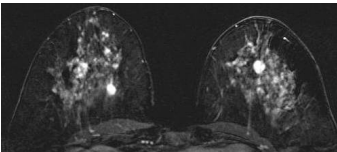

Subtraction image

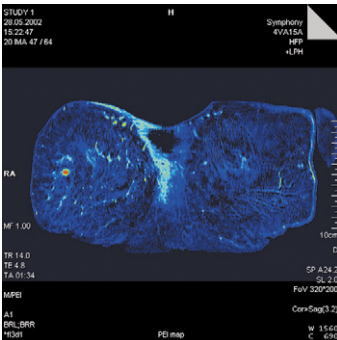


PEI parameter image

Performing motion correction (if required)

Measurement was performed MR images that are either not suitable for a diagnosis or only to a limited degree due to motion, can be corrected retroactively.

Two different motion correction algorithms are available for this purpose:

- “Fast” mode
Suitable for small to medium strong motion artifacts, faster computation times.
- “High Quality” mode
For very strong motion artifacts, however, with longer computation times.

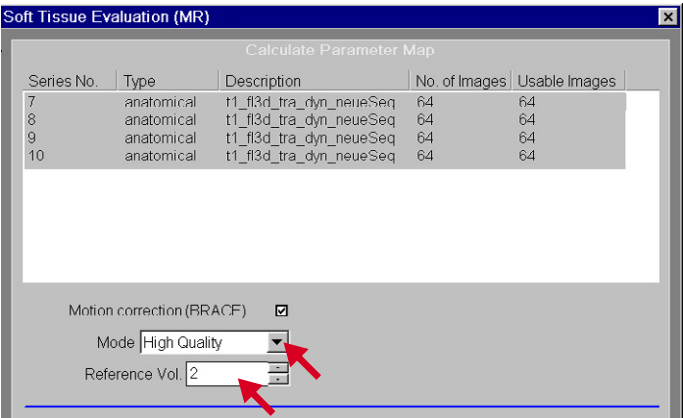

Perform motion correction



Select the desired motion correction mode.

Determine which measurement you want to use as a reference measurement for motion correction. (Reference volume)

TIP

The reference measurement should be the measurement where highest signal intensity can be expected in strongly enhanced tumors. This is the first post-contrast measurement for dynamic measurements with a temporal resolution between 60 to 90 seconds. For this reason select No. 2 as the reference volume.

Reconstructing combination images
(offline)

**Parameter
images are
available**

Optimum selection of all parameters is required to obtain a meaningful combination image. For this reason, only subsequent reconstruction of combination images is available.

Select e.g., a Washin and Washout series from the Patient Browser.

In the menu, select “Applications > Breast (MR)”.

Use the weighting factor to determine the extent to which each of the parameter images should contribute to the result.

The default value for Washout is – 1.0 and for Washin it is +1.0.

Effect of the weighting factor

Weighting factor	Parameter image	Combination image
negative	small pixel value	bright pixel
	large pixel value	dark pixel
positive	small pixel value	dark pixel
	large pixel value	bright pixel


Reconstruct combination images

Use the “Window” and “Center” parameters to set how the brightness and contrast of the parameter images will contribute to the combination image.

Please note: The following rule applies when using fixed values.

Pixel value $\leq$ Window – Center / 2 = lowest gray scale or color value
 Pixel value $\geq$ Window + Center / 2 = largest gray scale or color value

Start the reconstruction of combination images with “Evaluate”.

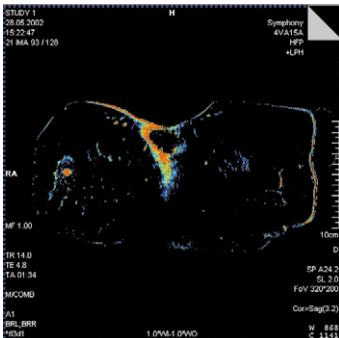

TIP

Take the optimum setting for brightness and contrast from the displayed windowing values in the Viewer (W/C) or set fixed values.

*Washin/
Washout
combination
image*

Evaluate contrast agent dynamics (offline)

Precontrast/ postcontrast images are available

Use the Mean Curve to evaluate the contrast agent dynamics.

Select all precontrast and postcontrast series in the Patient Browser.

In the menu select

"Applications > Mean Curve".

The images are loaded in the upper left segment.

In the Patient Browser, select the subtraction series in which you want to draw the ROI, and use drag & drop to move the image to the lower left segment. For improved localization of suspicious areas, you can also load an additional reference image (e.g., sagittal MIP image or parameter image) in the lower right segment.

Select the subtraction image you want to evaluate.

Copy this slice position to the upper left segment.

Draw the desired ROIs into the subtraction image.


The drawn ROIs are automatically transferred into the upper left segment.



Start the evaluation.

Confirm with "OK" that the pre-contrast series was taken into account as a reference series for the evaluation.

For results, you receive a curve evaluation and, behind it, a tabular evaluation, in the upper right segment.

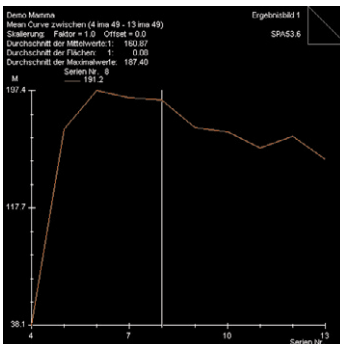

TIP

Storing the evaluation, including the anatomical image with the ROI:



Curve evaluation

Manual post-processing (offline)

If you have not selected inline reconstruction for the dynamic measurement, upon completion of the measurement you can reconstruct subtraction images from the series obtained and statistically analyze parameter images.

Manual subtraction

Select the precontrast and one or more postcontrast series in the Patient Browser.

In the menu, select “Evaluation > Dynamic Analysis > Subtraction”.

The precontrast series is automatically subtracted from the postcontrast series.

Click Test to compute an example subtraction image using the parameters currently set.

Click OK to generate the entire subtraction series.

The subtraction series is automatically stored in the Patient Browser.

Statistically analyze Washin and Washout images

In the Viewer, draw a ROI in the Washin or Washout image.

The minimum and maximum value in the ROI is displayed.

Divide the displayed values by 10. You obtain the minimum and maximum signal increase in the ROI in percent per minute.

TIP


If the minimum and maximum values are not displayed in the ROI, click in the ROI and open the context menu with the right mouse button.

Select “Minimum” and “Maximum” under menu item “Options”.

Manually scale subtraction images (offline)

Adaptation of the intensity range necessary

In some cases during subtraction, the result image may display “black” areas. This means the postcontrast data set has a lower signal intensity than the native data set. In this case, you can scale the intensity range for the subtraction.

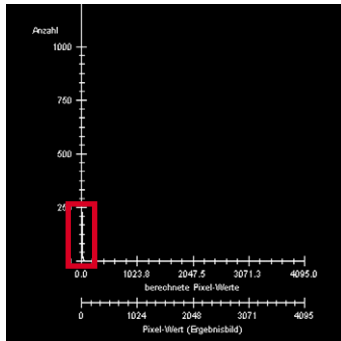


Scaling required:

The intensity values of the pixel are accumulated along the Y axis; i.e., they cannot be displayed in the intensity range.

Change the intensity range in the “Scaling” dialog box.

In this case, click the right mouse button in the histogram and then with the left mouse button click “Properties”.

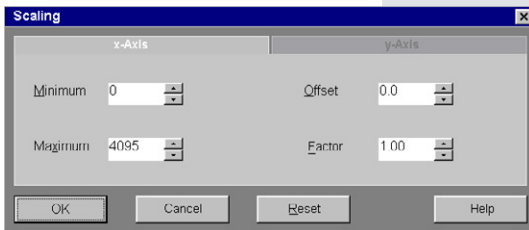

Scaling dialog box

Move the “X-Axis” function card to the foreground.


**Reset scaling,
e.g., Minimum –64,
Maximum 4032.**

The histogram is adjusted accordingly.

Prosthetic silicone implants

MR imaging of prosthetic silicone implants require images with high resolution. Depending on the clinical question, either fat or silicon can be suppressed, or pure silicon images can be created.

- | | |
|--------------------------|---|
| Measurement types | <ul style="list-style-type: none">• Pure silicon images:
to display capsule contractures, prosthesis dislocation, ruptures (intracapsular: “linguini” signs, extracapsular: defect in the fibrous capsule, possibly siliconoma).• Silicone suppression:
to display fat or tumors.• Water suppression:
to suppress the signal from cysts (water is displayed dark)• Fat suppression:
to display water-filled cysts (water is displayed bright). |
|--------------------------|---|

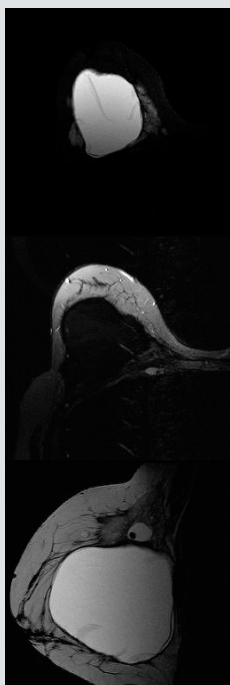


*pure
silicone image*

*with silicone
suppression*

*with water
presaturation*

Spectral differences

Spectral differences between fat, silicone, and normal breast parenchyma enable selective signal suppression.

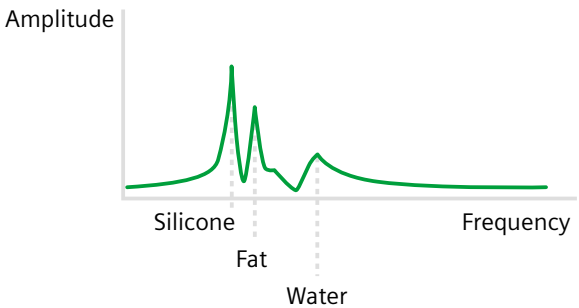

Display of MR spectra: The frequency of spectra (silicone, fat, and water) runs from left to right. Normally, the fat line somewhat overlays the silicone line.


Fat/silicone suppression

Fat or silicone suppression may be complete or partial.

Better orientation. Contours of the breast and implant may be better viewed using partial fat/silicone suppression.

Displaying pure silicone images

It may be necessary to display silicone implants during breast examinations. In this case, the signals from fat and water have to be suppressed. The result image displays the pure silicone signal.

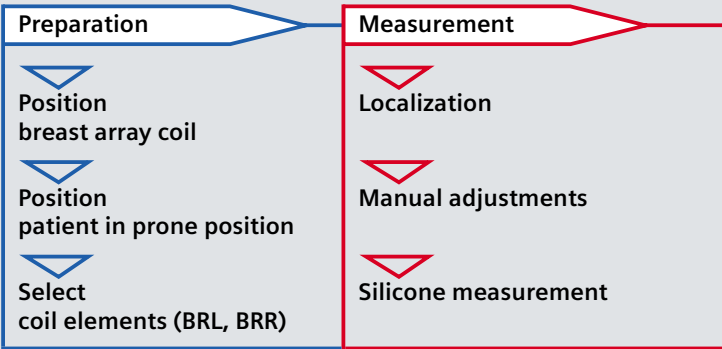

Recommendation. To obtain better adjustment results, we recommend these measurements be performed selectively on one breast only. As indicated under “Preparation”, select coil element BRL or BRR an (depending on the breast in question).

Possible measurements

- STIR with water saturation :
the fat signal is suppressed by an additional inversion pulse. The water saturation pulse (centered on the water peak) suppresses the signal from the breast and blood vessels. In the result image, silicone is displayed bright.
- T2 sagittal :
to check the impermeability of the silicone implant. Displays creases in the implant (e.g., linguini or keyhole signs)
- Silicone suppression :
with any T1-weighted measurement or T2-weighted measurement using a water saturation pulse (centered on the silicone peak). In the result image, silicone is displayed dark.

Spectral water saturation is possible because the frequency separation between water and silicone is greater than that between water and fat. However, silicone might become partially saturated.

Perform manual adjustments (I)

Protocol is selected


Position the slices for the respective protocol (fat saturation, water suppression, etc.)

On the System Adjustments parameter card, select “Confirm freq. adjustment”.

Start the measurement.

Adjustment is performed automatically.
The “Confirm Frequency Spectrum” dialog window is displayed.

*Peaks from left to right:
silicone—
fat—water*

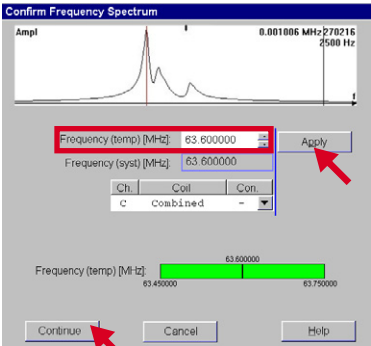

Click “Apply” and “Continue”, to copy the desired frequency.

If you want to change the frequency:

Enlarge the frequency spectrum (double-click with the *right* mouse button).

Improving fat saturation

Center the frequency on the fat peak (click the left mouse button on the fat peak).

This field is then copied to the field "Frequency (temp)".

Add 220 Hz at 1.5 Tesla to this value or 440 Hz at 3.0 Tesla (the last 3 numbers of this frequency).

Click "Apply" and "Continue".

The measurement is performed without additional adjustments.

TIP

With very high fat proportions: Check the automatic peak detection prior to each measurement with spectral fat sat and adjust manually, if necessary.

Perform manual adjustments (II)

Protocol is selected

Change water saturation


The water saturation suppresses the peak used for centering, and may be utilized for almost every protocol.

On the contrast card of the selected protocol, select “Water Suppression”.

Select the system card for the protocol.

Select parameters

“Confirm freq. adjustment” and “Assume Silicone”.

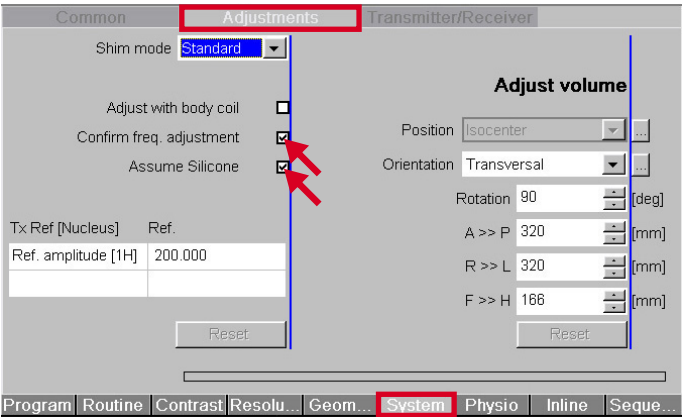

As soon as the “Confirm Frequency Spectrum” dialog window is displayed:

Click “Apply” and “Continue”, to copy the desired frequency.

If you want to change the frequency:

Click the left mouse button on the frequency peak you want to suppress.

TIP

For unilateral or small implants: Check the automatic peak detection, because the silicone signal can be very small. For this purpose, view the signals of the individual coil elements so you can check whether the system frequency on the side of the implant is located on the desired peak (e.g., water, if water is to be saturated to obtain a silicone image).

A	ADC map	22
	ADC value	22
B	Breast lesions	
	Combination images	11
	Dynamic examination	
	overview	18
	MIP images	7
	Motion correction	8, 30
	Parameter images	10
	Perform dynamic	
	measurement	28
	Prepare contrast agent	
	injection	20
	Prepare patient	20
	Qualitative evaluation	6, 7
	Quantitative evaluation	10, 11
	Reconstructing	
	combination images	32
	Select parameter images	26
	Sensitivity	4
	Set dynamic measurement	24
	Set subtraction	25
	Specificity	5
	Subtraction images	6
C	Combination images	11
	Reconstructing	32
	Weighting factor	32
	Contrast agent injection	
	Prepare	20

D	Diffusion-weighted measurement with REVEAL	22
	Dynamic examination	
	Mean Curve	34
	Overview	18
	Perform measurement	28
	Reconstructing combination images	32
	Select parameter images	26
	Set measurement	24
	Set subtraction	25
F	Fat suppression	41
H	Highest value (Washin)	12
M	Manual adjustment	46, 48
	Manual subtraction	36
	Mean Curve	34
	MIP images	7
	MIPT	10, 16
	Motion correction	8
	Performing it	30
	MR Mammography	
	Application	2
	Overview	2
P	Parameter images	10
	MIPT	10, 16
	PEI	10, 14
	Selecting	26
	TTP	10, 14
	Washin	10, 12
	Washout	10, 12

Patient	
Prepare for examination	20
PEI	10, 14
Prosthetic silicone implants	40
Displaying silicone images	44
Fat suppression	41
Manual adjustment	46, 48
Silicone images	41
Silicone suppression	41
Spectral differences	42
R REVEAL	22
S Scaling	
Manually scale	
subtraction images	38
Sensitivity	4
Silicon implants	40
Silicone images	41
Silicone suppression	41
Silicone images	41
Displaying	44
Silicone implants	
Displaying silicone images	44
Fat suppression	41
Manual adjustment	46, 48
Spectral differences	42
Silicone suppression	41
Specificity	5
Spectral differences	42
Statistical analysis	
Washin/Washout	37

	Subtraction images	6
	Manual scaling	38
	Manual subtraction	36
	Set subtraction	25
T	TTP	10, 14
W	Washin	10, 12
	Highest value	12
	Statistically analyze	37
	Washout	10, 12
	Statistically analyze	37
	Weighting factor	
	(combination images)	32

© Siemens AG 2008–2009
Order number
MR-07002.643.09.01.02
Printed in Germany
05/2009

Siemens AG
Wittelsbacherplatz 2
D-80333 Muenchen
Germany

Contact:
Siemens AG
Healthcare Sector
Henkestrasse 127
D-91052 Erlangen
Germany
Telefon: +49 9131 84-0
www.siemens.com/healthcare

www.siemens.com/healthcare