

GE Healthcare

LightSpeed 5.X RT

Pre-Installation Manual



OPERATING DOCUMENTATION



2384867-100
Rev 8

This Page Intentionally Blank.

IMPORTANT PRECAUTIONS

LANGUAGE

警告 (ZH-CN)	• 本维修手册仅提供英文版本。 • 如果维修服务提供商需要非英文版本，客户需自行提供翻译服务。 • 未详细阅读和完全理解本维修手册之前，不得进行维修。 • 忽略本警告可能对维修人员，操作员或患者造成触电、机械伤害或其他形式的伤害。
ПРЕДУПРЕЖДЕНИЕ (BG)	• Това упътване за работа е налично само на английски език. • Ако доставчикът на услугата на клиента изиска друг език, задължение на клиента е да осигури превод. • Не използвайте оборудването, преди да сте се консултирали и разбрали упътването за работа. • Неспазването на това предупреждение може да доведе до нараняване на доставчика на услугата, оператора или пациент в резултат на токов удар или механична или друга опасност.
VÝSTRAHA (CS)	• Tento provozní návod existuje pouze v anglickém jazyce. • V případě, že externí služba zákazníkům potřebuje návod v jiném jazyce, je zajištění překladu do odpovídajícího jazyka úkolem zákazníka. • Nesnažte se o údržbu tohoto zařízení, aniž byste si přečetli tento provozní návod a pochopili jeho obsah. • V případě nedodržování této výstrahy může dojít k poranění pracovníka prodejního servisu, obslužného personálu nebo pacientů vlivem elektrického proudu, respektive vlivem mechanických či jiných rizik.
ADVARSEL (DA)	• Denne servicemanual findes kun på engelsk. • Hvis en kundes tekniker har brug for et andet sprog end engelsk, er det kundens ansvar at sørge for oversættelse. • Forsøg ikke at servicere udstyret medmindre denne servicemanual har været konsulteret og er forstået. • Manglende overholdelse af denne advarsel kan medføre skade på grund af elektrisk, mekanisk eller anden fare for teknikeren, operatøren eller patienten.
WAARSCHUWING (NL)	• Deze onderhoudshandleiding is enkel in het Engels verkrijgbaar. • Als het onderhoudspersoneel een andere taal vereist, dan is de klant verantwoordelijk voor de vertaling ervan. • Probeer de apparatuur niet te onderhouden voordat deze onderhoudshandleiding werd geraadpleegd en begrepen is. • Indien deze waarschuwing niet wordt opgevolgd, zou het onderhoudspersoneel, de operator of een patiënt gewond kunnen raken als gevolg van een elektrische schok, mechanische of andere gevaren.

WARNING (EN)	• This Service Manual is available in English only. • If a customer's service provider requires a language other than English, it is the customer's responsibility to provide translation services. • Do not attempt to service the equipment unless this service manual has been consulted and is understood. • Failure to heed this warning may result in injury to the service provider, operator, or patient, from electric shock or from mechanical or other hazards.
HOIATUS (ET)	• Käesolev teenindusjuhend on saadaval ainult inglise keeles. • Kui klienditeeninduse osutaja nõuab juhendit inglise keelest erinevas keeles, vastutab klient tõlketeenuse osutamise eest. • Ärge üritage seadmeid teenindada enne eelnevalt käesoleva teenindusjuhendiga tutvumist ja sellest aru saamist. • Käesoleva hoiatuse eiramine võib põhjustada teenuseosutaja, operaatori või patsiendi vigastamist elektrilöögi, mehaanilise või muu ohu tagajärjel.
VAROITUS (FI)	• Tämä huolto-ohje on saatavilla vain englanniksi. • Jos asiakkaan huoltohenkilöstö vaatii muuta kuin englanninkielistä materiaalia, tarvittavan käännöksen hankkiminen on asiakkaan vastuulla. • Älä yritä korjata laitteistoa ennen kuin olet varmasti lukenut ja ymmärtänyt tämän huolto-ohjeen. • Mikäli tätä varoitusta ei noudateta, seurauksena voi olla huoltohenkilöstön, laitteiston käyttäjän tai potilaan vahingoittuminen sähköiskun, mekaanisen vian tai muun vaaratilanteen vuoksi.
ATTENTION (FR)	• Ce manuel de service n'est disponible qu'en anglais. • Si le technicien du client a besoin de ce manuel dans une autre langue que l'anglais, c'est au client qu'il incombe de le faire traduire. • Ne pas tenter d'intervenir sur les équipements tant que le manuel service n'a pas été consulté et compris. • Le non-respect de cet avertissement peut entraîner chez le technicien, l'opérateur ou le patient des blessures dues à des dangers électriques, mécaniques ou autres.
WARNUNG (DE)	• Diese Serviceanleitung existiert nur in Englischer Sprache. • Falls ein fremder Kundendienst eine andere Sprache benötigt, ist es aufgabe des Kunden für eine Entsprechende Übersetzung zu sorgen. • Versuchen Sie nicht diese Anlage zu warten, ohne diese Serviceanleitung gelesen und verstanden zu haben. • Wird diese Warnung nicht beachtet, so kann es zu Verletzungen des Kundendiensttechnikers, des Bedieners oder des Patienten durch Stromschläge, Mechanische oder Sonstige gefahren kommen.
ΠΡΟΕΙΔΟΠΟΙΗΣΗ (EL)	• Το παρόν εγχειρίδιο σέρβις διατίθεται στα αγγλικά μόνο. • Εάν το άτομο παροχής σέρβις ενός πελάτη απαιτεί το παρόν εγχειρίδιο σε γλώσσα εκτός των αγγλικών, αποτελεί ευθύνη του πελάτη να παρέχει υπηρεσίες μετάφρασης. • Μην επιχειρήσετε την εκτέλεση εργασιών σέρβις στον εξοπλισμό εκτός εάν έχετε συμβουλευτεί και έχετε κατανοήσει το παρόν εγχειρίδιο σέρβις. • Εάν δε λάβετε υπόψη την προειδοποίηση αυτή, ενδέχεται να προκληθεί τραυματισμός στο άτομο παροχής σέρβις, στο χειριστή ή στον ασθενή από ηλεκτροπληξία, μηχανικούς ή άλλους κινδύνους.

FIGYELMEZTETÉS (HU)	• Ezen karbantartási kézikönyv kizárólag angol nyelven érhető el. • Ha a vevő szolgáltatója angoltól eltérő nyelvre tart igényt, akkor a vevő felelőssége a fordítás elkészítése. • Ne próbálja elkezdni használni a berendezést, amíg a karbantartási kézikönyvben leírtakat nem értelmezték. • Ezen figyelmeztetés figyelmen kívül hagyása a szolgáltató, működtető vagy a beteg áramütés, mechanikai vagy egyéb veszélyhelyzet miatti sérülését eredményezheti.
AÐVÖRUN (IS)	• Þessi þjónustuhandbók er eingöngu fáanleg á ensku. • Ef að þjónustuveitandi viðskiptamanns þarfnast annas tungumáls en ensku, er það skylda viðskiptamanns að skaffa tungumálþjónustu. • Reynið ekki að afgreiða tækið nema að þessi þjónustuhandbók hefur verið skoðuð og skilin. • Brot á sinna þessari aðvörun getur leitt til meiðsla á þjónustuveitanda, stjórnanda eða sjúklings frá raflosti, vélrænu eða öðrum áhættum.
AVVERTENZA (IT)	• Il presente manuale di manutenzione è disponibile soltanto in inglese. • Se un addetto alla manutenzione richiede il manuale in una lingua diversa, il cliente è tenuto a provvedere direttamente alla traduzione. • Si proceda alla manutenzione dell'apparecchiatura solo dopo aver consultato il presente manuale ed averne compreso il contenuto • Il non rispetto della presente avvertenza potrebbe far compiere operazioni da cui derivino lesioni all'addetto, alla manutenzione, all'utilizzatore ed al paziente per folgorazione elettrica, per urti meccanici od altri rischi.
警告 (JA)	• このサービスマニュアルには英語版しかありません。 • サービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。 • このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないでください。 • この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。
경고 (KO)	• 본 서비스 지침서는 영어로만 이용하실 수 있습니다 . • 고객의 서비스 제공자가 영어 이외의 언어를 요구할 경우 , 번역 서비스를 제공하는 것은 고객의 책임입니다 . • 본 서비스 지침서를 참고했고 이해하지 않는 한은 해당 장비를 수리하려고 시도하지 마십시오 . • 이 경고에 유의하지 않으면 전기 쇼크 , 기계상의 혹은 다른 위험으로부터 서비스 제공자 , 운영자 혹은 환자에게 위해를 가할 수 있습니다 .
BRDINJUMS (LV)	• Šī apkāpes rokasgrāmata ir pieejama tikai angļu valodā. • Ja klienta apkāpes sniedzējam nepieciešama informācija citā valodā, nevis angļu, klienta pienākums ir nodrošināt tulkošanu. • Neveiciet aprīkojuma apkāpi bez apkāpes rokasgrāmatas izlasīšanas un saprašanas. • Šī brīdinājuma neievērošana var radīt elektriskās strāvas trieciena, mehānisku vai citu risku izraisītu traumā apkāpes sniedzējam, operatoram vai pacientam.

ĮSPĖJIMAS (LT)	• Šis eksploatavimo vadovas yra prieinamas tik anglų kalba. • Jei kliento paslaugų tiekėjas reikalauja vadovo kita kalba – ne anglų, numatyti vertimo paslaugas yra kliento atsakomybė. • Nemėginkite atlikti įrangos techninės priežiūros, nebent atsižvelgėte į šį eksploatavimo vadovą ir jį supratote. • Jei neatkreipsite dėmesio į šį perspėjimą, galimi sužalojimai dėl elektros šoko, mechaninių ar kitų pavojų paslaugų tiekėjui, operatoriui ar pacientui.
ADVARSEL (NO)	• Denne servicehåndboken finnes bare på engelsk. • Hvis kundens serviceleverandør trenger et annet språk, er det kundens ansvar å sørge for oversettelse. • Ikke forsøk å reparere utstyret uten at denne servicehåndboken er lest og forstått. • Manglende hensyn til denne advarselen kan føre til at serviceleverandøren, operatøren eller pasienten skades på grunn av elektrisk støt, mekaniske eller andre farer.
OSTRZEŻENIE (PL)	• Niniejszy podręcznik serwisowy dostępny jest jedynie w języku angielskim. • Jeśli dostawca usług klienta wymaga języka innego niż angielski, zapewnienie usługi tłumaczenia jest obowiązkiem klienta. • Nie próbować serwisować wyposażenia bez zapoznania się i zrozumienia niniejszego podręcznika serwisowego. • Niezastosowanie się do tego ostrzeżenia może spowodować urazy dostawcy usług, operatora lub pacjenta w wyniku porażenia elektrycznego, zagrożenia mechanicznego bądź innego.
ATENÇÃO (PT)	• Este manual de assistência técnica só se encontra disponível em inglês. • Se qualquer outro serviço de assistência técnica solicitar estes manuais noutro idioma, é da responsabilidade do cliente fornecer os serviços de tradução. • Não tente consertar o equipamento sem ter consultado e compreendido este manual de assistência técnica. • O não cumprimento deste aviso pode pôr em perigo a segurança do técnico, do operador ou do paciente devido a choques elétricos, mecânicos ou outros.
ATENȚIE (RO)	• Acest manual de service este disponibil numai în limba engleză. • Dacă un furnizor de servicii pentru clienți necesită o altă limbă decât cea engleză, este de datoria clientului să furnizeze o traducere. • Nu încercați să reparați echipamentul decât ulterior consultării și înțelegerii acestui manual de service. • Ignorarea acestui avertisment ar putea duce la rănirea depanatorului, operatorului sau pacientului în urma pericolelor de electrocutare, mecanice sau de altă natură.
ОСТОРОЖНО! (RU)	• Данное руководство по обслуживанию предлагается только на английском языке. • Если сервисному персоналу клиента необходимо руководство не на английском, а на каком-то другом языке, клиенту следует самостоятельно обеспечить перевод. • Перед обслуживанием оборудования обязательно обратитесь к данному руководству и поймите изложенные в нем сведения. • Несоблюдение требований данного предупреждения может привести к тому, что специалист по обслуживанию, оператор или пациент получат удар электрическим током, механическую травму или другое повреждение.

UPOZORNENIE (SK)	• Tento návod na obsluhu je k dispozícii len v angličtine. • Ak zákazníkovi poskytovateľ služieb vyžaduje iný jazyk ako angličtinu, poskytnutie prekladateľských služieb je zodpovednosťou zákazníka. • Nepokúšajte sa o obsluhu zariadenia skôr, ako si neprečítate návod na obsluhu a neporozumiete mu. • Zanedbanie tohto upozornenia môže vyústiť do zranenia poskytovateľa služieb, obsluhujúcej osoby alebo pacienta elektrickým prúdom, do mechanického alebo iného nebezpečenstva.
ATENCION (ES)	• Este manual de servicio sólo existe en inglés. • Si el encargado de mantenimiento de un cliente necesita un idioma que no sea el inglés, el cliente deberá encargarse de la traducción del manual. • No se deberá dar servicio técnico al equipo, sin haber consultado y comprendido este manual de servicio. • La no observancia del presente aviso puede dar lugar a que el proveedor de servicios, el operador o el paciente sufran lesiones provocadas por causas eléctricas, mecánicas o de otra naturaleza.
VARNING (SV)	• Den här servicehandboken finns bara tillgänglig på engelska. • Om en kunds servicetekniker har behov av ett annat språk än engelska ansvarar kunden för att tillhandahålla översättningstjänster. • Försök inte utföra service på utrustningen om du inte har läst och förstår den här servicehandboken. • Om du inte tar hänsyn till den här varningen kan det resultera i skador på serviceteknikern, operatören eller patienten till följd av elektriska stötar, mekaniska faror eller andra faror.
DIKKAT (TR)	• Bu servis kılavuzunun sadece ingilizcesi mevcuttur. • Eğer müşteri teknisyeni bu kılavuzu ingilizce dışında bir başka lisandan talep ederse, bunu tercüme ettirmek müşteriye düşer. • Servis kılavuzunu okuyup anlamadan ekipmanlara müdahale etmeyiniz. • Bu uyarıya uyulmaması, elektrik, mekanik veya diğer tehlikelerden dolayı teknisyen, operatör veya hastanın yaralanmasına yol açabilir.

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "damage in shipment" written on all copies of the freight or express bill before delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 8.
- Fill out a report on <http://egems.med.ge.com/edq/home.jsp>
- Contact your local service coordinator for more information on this process.

Rev. June 13, 2006

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Healthcare personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

IMPORTANT...X-RAY PROTECTION

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Healthcare, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protection, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Healthcare, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

Various protective materials and devices are available. It is urged that such materials or devices be used.

LITHIUM BATTERY CAUTIONARY STATEMENTS



CAUTION
Risk of
Explosion

Danger of explosion if battery is incorrectly replaced. Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer's instructions.



ATTENTION
Danger
d'Explosion

Il y a danger d'explosion s'il y a remplacement incorrect de la batterie. Remplacer uniquement avec une batterie du même type ou d'un type recommandé par le constructeur. Mettre au rebut les batteries usagées conformément aux instructions du fabricant.

OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the GEMS CQA Process to report all omissions, errors, and defects in this publication.

Revision History

Revision	Date	Reason for change
8	12/04/07	Chapter 4: - Section 1.0: Added bullet point concerning cable storage. - Section 2.1.1: Corrected Regulatory Minimum Working Clearance by Major Subsystem tables to match what is in the LS 5.X Installation Manual. Added tables with information for the UPS and A1 Disconnect. - Section 4.1: Added bullet point concerning clear access for service components. - Section 8.4 Floor levelness: added Note about specifications for the floor/gantry to prevent cable crushing. - Section 9.0: Added the network specifications. Chapter 5: Section 4.0: Added note concerning EMI testing near power sub-stations. Chapter 9: Section 7.1: Added information concerning LOTO and UPS.
7	08/22/07	Added Bulgarian Warning in Important Precautions Chapter 4 and Chapter 6 Updated Gantry Dimension Chapter 4: Added Limited Access contents

Revision	Date	Reason for change
6	04/19/07	Chapter 1: Added Section 4 - Medical Electrical Equipment for EMC. Chapter 2: - Moved Checklist and Block Diagram to Chapter 1. - Added Notice to review DVD prior to beginning Pre-Install work. - Added "Floor Levelness Specification" Chapter 4: Major updates, including (but not limited to) definition of <i>Regulatory</i> vs. <i>Service</i> Clearances and updates to "Structural Requirements". Chapter 5: - Updated Section 1.0 - Temperature and Humidity Specifications - Updated Section 2.0 - Cooling Requirement - Delete Section 4.9 EMC Edition 2 Compliance - duplicate with newly added Chapter 1, Section 4. Chapter 6: - Added Table 6-2 - Gantry and Table Mounting Requirements - Updated Figure 6-5 - Power Distribution Unit (NGPDU) Chapter 7: - Added Notice regarding lifting of CT Gantry. - Updated Section 4.0 - Storage Requirements - Updated Section 5.0 - Extreme Temperature Transportation and Deliveries - Updated Section 6.0 - System Transportation Chapter 8: - Updated Section 2.1 Power Source Configuration for PQR 13106373 Chapter 9: - Updated Figure 9-1 System Interconnect Diagram for include options. - Added Section 3.3 cable list of options - Added Section 5 Fuse information - Added Notes in Section 1.0: Use dry cleaning for electro components. - Updated Section 4.0 - Contractor Supplied Components - Updated Section 6.0 - UPS Interconnect - Updated Section 7.0 - Typical Customer Supplied Wiring Delete Appendix "FCT Installation Site Ready Form" Updated formatting to latest "branding" standards.
5	08/28/05	PSR 13050339: Added Note to Figure 9-2: Interconnection Runs
4	07/21/05	PSR 13033111 - Changed Scatter Survey values (shown in Chapter 4) from mR/scan to µGy/scan.
3	10/22/04	EMC Ed2 Update. Chapter 6 - Updated Floor Anchor information. Chapter 9 - Added UPS Compatability Matrix.

Revision	Date	Reason for change
2	05/18/04	CTCge87943: Corrected Floor Anchoring Information. Chapter 5 - Updated Cooling Requirements Worksheet. Chapter 6 - Corrected Gantry Floor Loading Spec.
1	12/15/03	Chapter 5 - Updated Cooling Requirements
0	09-15-03	Initial Release

Table of Figures

Figure 2-1	Hole Locations	43
Figure 3-1	Base Scanner System.....	46
Figure 4-1	Regulatory Clearance Requirements for LightSpeed RT System Configuration	48
Figure 4-2	Room Layout with cabinets	59
Figure 4-3	Room Layout without cabinets	59
Figure 4-4	Minimum Service Clearance	60
Figure 4-5	Typical Control Room Layout.....	63
Figure 4-6	Typical Room Layout with CT Simulation Laser Alignment System	64
Figure 4-7	Gantry Front Cover with Service Dolly Dimensions	65
Figure 4-8	Table and Gantry (Side View)	66
Figure 4-9	Gantry shown tilted +30° (top) and -30° (bottom).....	67
Figure 4-10	Power Distribution Unit (NGPDU-3)	68
Figure 4-11	Power Distribution Unit (NGPDU-4)	68
Figure 4-12	Operator's Console	69
Figure 4-13	Typical Room Layout, with Injector Option.....	71
Figure 4-14	Console Rear Bulkhead	75
Figure 4-15	Typical Scatter Survey (Body Phantom)	79
Figure 4-16	Typical Scatter Survey (Head Phantom).....	80
Figure 5-1	Sample Room Layout, showing approximate EMI requirements	85
Figure 6-1	Typical Floor Anchor, Gantry and Table	89
Figure 6-2	Gantry Center of Gravity	90
Figure 6-3	Gantry Anchor Locations.....	90
Figure 6-4	Patient Table (CT1).....	91
Figure 6-5	Power Distribution Unit (NGPDU)	92
Figure 6-6	NGPDU Seismic Mounting Bracket (dimensions are in mm)	93
Figure 6-7	Operator's Console	94

Figure 7-1	
Gantry Strap Location.....	95
Figure 7-2	
Gantry with Shipping Dollies and Side Rails.....	99
Figure 7-3	
Console ready to be unloaded from shipping skid.....	100
Figure 8-1	
System Ground Map.....	105
Figure 9-1	
System Interconnect Diagram	109
Figure 9-2	
Interconnection Runs.....	113
Figure 9-3	
Typical UPS Interconnect	116
Figure 9-4	
Primary Power Disconnect (A1) - Typical GE A1 shown.....	117
Figure 9-5	
Primary Power Disconnect (A1) – Fusible Disconnect and Magnetic Contactor	118
Figure 9-6	
TS6 X-Ray Warning Light Connections	119
Figure 9-7	
TS6 Room Door Interlock Connections - Without a Door Interlock	119
Figure 9-8	
TS6 Room Door Interlock Connections - With a Door Interlock	119

Table of Tables

Table 1-1	Schedule Date Commitments	31
Table 1-2	General Site Planning	32
Table 1-3	References for Specific Site Requirements	32
Table 1-4	Equipment Compatibility	33
Table 1-5	Network Connections	33
Table 1-6	Miscellaneous Tasks	33
Table 1-7	Emission Declaration	35
Table 1-8	Immunity Declaration	36
Table 1-9	Immunity Declaration con't.....	37
Table 1-10	Separation Distances	39
Table 3-1	System Options Catalog (Part) Numbers	45
Table 4-1	Console Subsystem	49
Table 4-2	NGPDU Subsystem	49
Table 4-3	Gantry Subsystem.....	50
Table 4-4	Table Subsystem	50
Table 4-5	UPS Subsystem	51
Table 4-6	A1 Disconnect Subsystem	52
Table 4-7	Equipment to be stored in storage cabinet.....	62
Table 4-8	Dimensions of Components	65
Table 4-9	Minimum Dimensions and Clearances.....	70
Table 4-10	Shielding Requirements Scaling	78
Table 5-1	Cooling Requirements (Worksheet)	82
Table 6-1	LightSpeed RT System Floor Loading	87
Table 6-2	Gantry and Table Mounting Requirements	88
Table 7-1	Size of Gantry & Dollies, with and without Side Rails	99

Table 8-1

Nominal Line Voltage

103

Table 8-2

Purchasable Options

103

Table 8-3

Minimum Feeder Wire Size

103

Table 8-4

Minimum Sub-Feeder Wire Size.....

103

Table 9-1

Component Designators

107

Table 9-2

GE Healthcare Supplied Cables (Standard Run) - UL Information.....

108

Table 9-3

GE Healthcare Supplied Cables (Optional, Long Run) - UL Information.....

110

Table 9-4

GE Healthcare Supplied Cables for Options - UL Information.....

111

Table 9-5

Runs 1, 2, 3, 4 and 5 Connections

112

Table 9-6

Contractor-Supplied Components

114

Table 9-7

H-Power Fuse Kit (2218570-2 BOM, rev 2).....

115

Table 9-8

LightSpeed RT Partial UPS Back-up Options.....

116

Table of Contents

Chapter 1

Introduction..... 25

Section 1.0

Site Readiness..... 25

Section 2.0

Responsibility of Purchaser..... 26

2.1 Customer Room Prep Items 26

2.2 Purchaser Site Preparation Work 27

2.3 Manufacturer's System Level Siting Requirements 27

2.3.1 Meeting Site Ready Requirements 28

2.3.2 Quick Installs 28

2.3.3 "Two-Step" and Upgrade Installs 28

2.3.4 Site Ready Inspection Visit 29

Section 3.0

Pre-Installation Checklist 31

Section 4.0

Medical Electrical Equipment for EMC..... 35

4.1 General Scope 35

4.2 Electromagnetic Emission 35

4.3 Electromagnetic Immunity 36

4.3.1 Limitations Management : 39

4.4 Use Limitation : 39

4.4.1 External components 39

4.5 Installation Requirements & Environment Control : 40

4.5.1 Cable Shielding & Grounding 40

4.5.2 This product complies the radiated emission as per CISPR11 Group 1 Class A standard limits 40

4.5.3 Subsystem & Accessories Power Supply Distribution 40

4.5.4 Stacked Components & Equipment..... 40

4.5.5 Low Frequency Magnetic Field..... 40

4.5.6 Static Magnetic Field Limits 40

4.5.7 Electrostatic Discharge Environment & Recommendations 40

Chapter 2

Pre-Installation Overview..... 41

Section 1.0

Dust/Dirt Contamination 41

Section 2.0

Chemical Contamination 41

Section 3.0

Walls, Ceiling, and Floor..... 41

Section 4.0	
Broad-band.....	42
Section 5.0	
Phone Line (for optional modem)	42
Section 6.0	
Review	42
 Chapter 3	
System Catalog	45
Section 1.0	
Option Catalog Numbers	45
Section 2.0	
Base Scanner System	46
2.1 Application.....	46
2.2 Configuration	46
 Chapter 4	
Room Planning	47
Section 1.0	
Required Systems Clearances	47
Section 2.0	
Regulatory and Service Clearances.....	48
2.1 Regulatory Clearances.....	48
2.1.1 Regulated Minimum Working Clearance by Major Subsystem	49
2.1.2 Terms and Definitions	53
Section 3.0	
Additional Regulatory Clearance Information	54
3.1 Minimum Room Size (Limited Access)	54
3.1.1 Regulatory Caution	54
3.1.2 Egress Clearance.....	54
3.1.3 Operational Caution	55
3.1.4 System Specifications (LightSpeed RT Systems).....	55
3.1.4.1 Recommended Room Size	55
3.1.4.2 Typical Room Size	55
3.1.4.3 Minimum Room Size.....	55
3.2 How to Measure	56
3.3 Minimum Room Size & Requirement Layouts	57
3.4 Recommended and Typical Room Size & Requirement Layouts	59
Section 4.0	
Service Clearance.....	60
4.1 Service Clearances for Single Service Engineer	61
4.2 Power Distribution Unit (NGPDU) Service Clearance.....	61
4.3 Console Service Clearance.....	61

Section 5.0	
Recommended Layouts	62
5.1 Storage Cabinet.....	62
5.2 Advantage Workstation (AW)	62
5.3 Control Room Considerations.....	63
5.4 CT Simulation Laser Alignment Systems	64
Section 6.0	
Component Dimensions	65
6.1 Table and Gantry	66
6.2 Power Distribution Unit	68
6.3 Operator's Console.....	69
Section 7.0	
Minimum Dimensions and Clearances.....	70
7.1 Injector Control	70
Section 8.0	
Structural Requirements	72
8.1 Table and Gantry Mounting Requirements.....	72
8.2 Floor Anchors	72
8.3 Floor Strength.....	73
8.4 Floor Levelness	73
8.5 Floor Vibration	73
8.5.1 Steady State Vibration.....	73
8.5.2 Transient Vibration.....	73
8.5.3 Equipment Location	73
8.6 Walls: Scan Window	74
Section 9.0	
Network Connections	74
9.1 US Broad Band Process Overview.....	75
9.2 Customer Broad-Band Responsibilities	75
Section 10.0	
Radiation Protection	78

Chapter 5	
Environmental Conditions	81
Section 1.0	
Temperature and Humidity Specifications	81
1.1 Temperature (Scan Room & Control Room)	81
1.2 Humidity (Scan Room & Control Room)	81
1.3 Other Guidelines.....	82
Section 2.0	
Cooling Requirement.....	82
Section 3.0	
Altitude.....	83

Section 4.0	
Electro-Magnetic Interference (EMI)	83
4.1 Gantry	83
4.2 Color Monitor.....	83
4.3 Console / Computer Equipment	83
4.4 Magnetic Media	83
4.5 PDU.....	83
4.6 EMI Reduction.....	83
4.7 UPS.....	84
4.8 Equipment EMI “Envelopes”	85
 Chapter 6	
Floor Loading and Weights	87
Section 1.0	
Floor Loading.....	87
Section 2.0	
Mounting Data, Including Seismic	88
 Chapter 7	
Delivery Data	95
Section 1.0	
Van Delivery	95
Section 2.0	
Delivery/Shipping Requirements & Considerations.....	95
Section 3.0	
Site Environmental Considerations	96
3.1 Dust/Dirt Contamination	96
3.2 Chemical Contamination	96
Section 4.0	
Storage Requirements	97
Section 5.0	
Extreme Temperature Transportation and Deliveries	97
Section 6.0	
System Transportation.....	98
Section 7.0	
Gantry Considerations	99
Section 8.0	
Table Considerations	100
Section 9.0	
Console Considerations	100

Chapter 8

Power Requirements 101

Section 1.0

Introduction 101

Section 2.0

System Input Power 101

2.1 Power Source Configuration 101

2.2 Rating 101

2.3 Regulation 102

2.4 Phase Imbalance 102

2.5 Sags, Surges & Transients 102

2.6 Grounding 102

Section 3.0

Recommended Power Distribution System 102

Section 4.0

Ground System 105

Chapter 9

Interconnection Data 107

Section 1.0

Introduction 107

Section 2.0

Component Designators 107

Section 3.0

Interconnect Runs, Wiring and Cables 108

3.1 GE Healthcare Supplied (Standard Length) 108

3.2 GE Healthcare Supplied (Optional, Long Run) 110

3.3 GE Healthcare Supplied (Cables of Options) 111

3.4 Contractor (Customer) Supplied 112

Section 4.0

Contractor Supplied Components 114

Section 5.0

Fuse 115

Section 6.0

UPS Interconnect 116

Section 7.0

Typical Customer Supplied Wiring 117

7.1 Primary Power Disconnect 117

7.2 Scan Room Warning Light & Door Interlock 119

End of Section.

Chapter 1

Introduction

This direction contains physical and electrical data necessary for planning and preparing a site. Pre-installation work is defined as site preparation for installation of the GE CT scanner. It is the responsibility of the purchaser to arrange and pay for this work. Pre-installation work includes:

- Installation of electrical conduit, junction boxes, ducts, outlets, and line safety switches.
- Installation of interconnection wiring that is AWG stranded copper. The electrical contractor shall ring out and tag all wires at both ends. Color-coded wires are recommended for easier identification. Wires shall be continuous without splices. Ground wires must conform to local codes.
- Any site renovation.
- Alterations and modifications to products not specifically included in the sales contract.

All work must conform to local building and safety codes. Unless specifically mentioned, GE Healthcare does not provide or install wires, conduits, junction boxes, and ducts as illustrated in this publication.

All CT site plans, preliminary concepts and final working drawings must be reviewed by General Electric Headquarters Architectural Planning prior to construction or approval.

Contact your local General Electric sales representative for complete information regarding your site-specific room layout.

Section 1.0

Site Readiness

Site ready is a requirement that must be achieved to install a CT product. For your convenience, a site ready visit inspection shall be performed at least three (3) days prior to the installation date. The site inspection must conclude with a minimum of a conditional pass status to be ready on the requested installation delivery date. Site ready inspections on the delivery date will not be acceptable unless prior arrangements have been made.

Pre-Installation and Site Ready Tools:

- Pre-Installation Video
- Floor template
- Pre-Installation check List
- Pre-Installation Block Diagram
- Site Room Layouts
- Power and Grounding Inspection
- Pre-Installation Support

Section 2.0

Responsibility of Purchaser

2.1 Customer Room Prep Items

The CT air intake is near the bottom of the gantry and draws air in through a filter in the gantry heater assembly. Fine dust as listed below will clog the filter and be deposited throughout the gantry, table, console and PDU electronics. This fine dust cannot be completely removed and can be damaging to electronic components.

For these reasons, the scanner should be the last item installed in your CT suite area.

“Pre-installation” is work necessary to plan and prepare a site for installation of equipment.

Pre-installation work helps the user (customer) avoid:

- Application delay and scheduling
- Surprise siting discoveries
- Installation confusion
- Waste of manpower

The following **MUST** be completed before installation work can begin for a GE CT scanner:

- Completely finished:
 - Wall painted or have final wall covering
 - Ceiling tiles installed and no remaining ceiling work is required
 - Final floors covering installed with no remaining dust causing floor work required
 - All room millwork installed as shown on the site print
 - All plumbing work in the CT suite is completed
 - No construction in or around the scan suite AREA that will produce:
 - * Concrete dust
 - * Drywall dust
 - * Ceiling tile dust
 - * Wood sawdust or shaving
 - * Dust tracked into the CT suite area
 - Active Broad Band connection
 - A completed network connection is required for ALL CT installations.
 - A GE Healthcare network specialist may be required to complete the VCN connection. This may take a week or longer to schedule.
 - Power available to A1, with provision for Lockout/Tagout at the A1 disconnect
- If a UPS is required, a GE A1 breaker* will be needed to complete this installation. Refer to the electrical section for more details.



NOTICE

SERVICE NOTICE: An improperly prepared site (i.e., one that is in a state of construction) can result in increased installation time.

A CT scanner installed in a dirty environment is more prone to contamination, which can result in decreased reliability and increased scanner downtime.

2.2 Purchaser Site Preparation Work

This list below will describe many of the items to consider when planning for a system replacement or designing a room for new equipment

- Determine room dimensions and verify that doorways are large enough for the scanner system.
- Install appropriate conduits and duct work for system cables. If additional components are required in the CT suite, their connection consideration must be determined and completed.
- Install junction boxes of correct size with covers at locations shown in installation plan.
- A1 main disconnect installation
- Install power supply of correct voltage output and adequate KVA rating.
- Install local disconnects, including proper over-current protection.
- Install “steelwork” or other suitable support work for mounting equipment on walls or from ceiling.
- Camera should be on-site at the time of installation.
- Complete all suite and room alterations and modifications.
- Verify that room shielding is adequate for the system being installed
- Review structural requirement - including floor vibration, levelness, and thickness
- Review HVAC requirements including system regulation and patient comfort.
- Review operational clearances to see if your daily used items fit, such as beds and carts.
- Emergency medical equipment should also be considered
- Storage cabinets and sink (if required) must be shown on the site print
- These contractors and others may be required to help confirm that the site meet all installation requirements:
 - Structural Engineer and /or Architect
 - HVAC contractor
 - Electrical contractor
 - Qualified radiological health physicist
- The project manager or sales will deliver a copy of the pre-installation quick start kit.

The above items can be found in chapter 2 through 9 in this manual.

It is suggested that this work be completed at least three days prior to delivery

2.3 Manufacturer's System Level Siting Requirements

These siting requirements are the minimum that must be met in order to install a new or replacement system.

- Network Communication in place and active
- Meets all scan room regulatory and service requirements
- Meets all minimum scan room structural requirements
- Meets minimum scan room HVAC requirements
- Meets minimum scan room electrical requirements
- Reviewed radiation protection section in the Pre-Installation manual
- All in room items shown on the final GE Healthcare site print and the final print is on site
- No construction in the scan room or neighboring suite areas

It is suggested that this work be completed at least three days prior to delivery

2.3.1 Meeting Site Ready Requirements

The site ready visit will take place at least three days prior to the delivery date. The site ready visit is intended to verify that all of the siting requirement are met and the site is ready for installation.

The site ready visit will result in a report to the project manager indicating one of the following:

Pass - All required items are present, completed and the site is ready for installation.

Conditional Pass - is issued when 80% of all of the task are completed and all parties agree that the 20% will be completed by the installation delivery date.

If a "Conditional Pass" is granted on the inspection date, the project manager must present conclusive evidence that unfinished tasks are completed and that the site is ready for delivery one business day prior to delivery.

Fail - is issued when less than 80% of the task are completed and all parties cannot agree that the remaining work will be completed by the requested installation delivery date. Failed sites will be rescheduled when all items are completed.

2.3.2 Quick Installs

Quick installations are described as sites with minimum room improvements required. These include, but are not limited to the following items:

- Existing electrical disconnect device, wire size and grounds meet all of the above requirements.
- Existing structural items including floor thickness meet all of the above requirements
- Existing HVAC capacity and regulation meet all of the above requirements
- Existing CT suite meets all of the above regulatory and minimum size requirements
- Existing facility can accommodate the delivery and meet all of the above delivery requirements

Quick Installs are subject to the following restrictions:

- Quick installs must have a new room print that accurately reflects the rooms to be upgraded.
- Existing floor anchors from a non-5.X system CANNOT be reused.
- New floor anchors must be a minimum of 4" from any existing floor penetrations.

Quick Installs typically involve a weekend de-install and room prep completion, with a next business day delivery and install.

2.3.3 "Two-Step" and Upgrade Installs

A "Two-Step" installation is the practice of temporarily installing one CT system in a site with the intention of upgrading the site to a different CT system at a later date.

- For a "two-step" installation to be considered, the room must meet the minimum room requirements for the project being upgraded.
- As with any upgrade installation, "two-steps" are subject to ALL of the siting requirements imposed by the upgrade/final system. This includes the recommended room size as well as electrical, structural and HVAC requirements.
- Two-steps and other upgrades may be done as "Quick" Installs. In this case, all requirements described in Section 2.3.2 (above) also apply.
- It is the customer's responsibility to check that all requirements are met.
- Rooms that do not meet the minimum requirements for the final product must either upgrade (or enlarge) their room, or consider the "Left-Side Limited Access" option.

2.3.4 Site Ready Inspection Visit

Must meet all in section 1 and 2 plus these additions reviewed at the site ready visit.

The GE Healthcare project manager will review the site delivery process with you to determine how to best transfer the equipment from the transportation truck to your room.

This site ready inspection will review and check these items:

Delivery information

- Determine delivery route into the scan room
- Determine if tilt dollies or riggers are required
- Determine if elevators, doorways and hallways are adequate for delivery
- Determine if floor protection is required
- Determine if a tilt bed truck is required for ground delivery and ordered.

Regulatory Requirement

- Room size meets the minimum requirements
- Site print is present and accurately reflects the room size and layout.
- No grounded walls are present in the regulatory clearance areas
- All regulatory clearances space is met
- Room meets all local codes

Manufacturer Requirements: As listed in section 2 - all requirements are met

Purchaser's Site Preparation Work: As listed in section 1 - all actions are completed

Section 3.0

Pre-Installation Checklist

Required Information for Site

Must be completed before the scheduled delivery date

Hospital Name as it appears on the system screens:

Network ID numbers / IP addresses

PACS: _____

Camera: _____

AW: _____

Other - Specify type & ID: _____

Other - Specify type & ID: _____

Camera setup information: _____

AW Direct Connect address: _____

Do you want HIPAA enabled? No___ Yes___

Do you want automatic downloads enabled? No___ Yes___

:

Table 1-1 Schedule Date Commitments

GE	Cust	Dates
Y N	Y N	
	☐ ☐	Has the project schedule been verified with facilities department, contractor, and GE?
	☐ ☐	Will the committed site-ready date be met?
☐ ☐	☐ ☐	Does the completion date for any/all construction meet or precede the delivery date?
☐ ☐	☐ ☐	Is the Power & Ground survey complete? Date: _____ Hospital contact: _____
	☐ ☐	Site-Ready visit is scheduled. Date: _____
	☐ ☐	Delivery date is scheduled. Date: _____
	☐ ☐	Installation date is scheduled. Date: _____
☐ ☐	☐ ☐	Installation timing: A: Weekdays___ B: Weekend___ C: Quick Install___ If B or C, have all sub-contractors been notified? No___ Yes___
☐ ☐	☐ ☐	Does the delivery and/or installation date need to be adjusted?
	☐ ☐	First-Use date is scheduled. Date: _____
☐ ☐	☐ ☐	Applications/Training dates: On-Site Training Date: _____ Healthcare Institute Training Date: _____

Table 1-2 General Site Planning

GE	CUST	General / Site Requirements
Y N	Y N	<i>Must be completed 5 weeks before scheduled delivery date</i>
☐ ☐	☐ ☐	Have final drawings been approved and distributed to the contractors?
☐ ☐	☐ ☐	Are final drawings “signed off” to approve equipment layout / orientation?
☐ ☐	☐ ☐	Do the actual room dimensions match those on the final drawings?
	☐ ☐	Has the radiologist health physician reviewed and approved the room layout and shielding requirements?
	☐ ☐	Have any additional requirements or questions about the installation been discussed with GE? List: _____ _____ _____ _____
	☐ ☐	Is there a person assigned to review and verify that all installation requirements are met? Name: _____
	☐ ☐	Have the specific site requirements been discussed with the contractor? Refer to the GE final drawings specifications. (See Table 1-3 below.)
	☐ ☐	Has the responsibility of cabling, installing, and interfacing accessories not on the order been discussed?
	☐ ☐	Are all third-party vendors identified, notified and scheduled? (Examples: Netcom, Medrad, etc.)
	☐ ☐	Have all regulatory requirements been met per <i>Regulatory and Service Clearances</i> , on page 48?
	☐ ☐	Will existing network, broadband, and camera cable drops reach new locations and will they meet the requirements and function with LightSpeed 5.X RT? If not, what are the requirements? List: _____ _____ _____

Table 1-3 References for Specific Site Requirements

Sections for Specific Requirements	
Room Planning - Chapter 4 (page 47)	Floor Loads & Weights - Chapter 6 (page 87)
Radiation Protection - page 78	Delivery - Chapter 7 (page 95)
Environment - Chapter 5 (page 81)	Power - Chapter 8 (page 101)
All work by contractors must be completed before the scheduled delivery date.	

Table 1-4 Equipment Compatibility

GE	Cust	Equipment
Y N	Y N	<i>Must be completed 5 weeks before scheduled delivery date</i>
☐ ☐	☐ ☐	Has the order been reviewed for completeness and compatibility with existing equipment? Typical equipment: Remote monitors ____ AW relocation ____ Cardiac option ____ Injectors ____
☐ ☐		Are interfaces to existing and/or new accessories ordered and planned for accordingly?
☐ ☐	☐ ☐	Have the following peripheral locations been included in the site drawings? EKG monitor ____ Injector control ____ Laser camera ____ UPS ____ 2 nd Monitor ____
☐ ☐	☐ ☐	Will GE Healthcare provide additional services per contract negotiations?
☐ ☐		Are correct length cables on order?

Table 1-5 Network Connections

GE	Cust	Network Installation and Setup
Y N	Y N	<i>Must be completed 5 weeks before scheduled delivery date</i>
☐ ☐	☐ ☐	Have IP addresses and Host Names been obtained? No ____ Yes ____
☐ ☐	☐ ☐	Will a network camera be used? No ____ Yes ____
☐ ☐	☐ ☐	Mandatory: Is the network installed, are network jacks installed, and is the entire network tested?
☐ ☐	☐ ☐	Mandatory: Broadband VPN installed/setup?
☐ ☐	☐ ☐	Mandatory: Are network software options ordered ____ HIS RIS option ____ DICOM print ____ AW ____
	☐ ☐	Optional: Has modem option ordered? ____ (Requires a site escalation)
☐ ☐	☐ ☐	Optional: Is the LightSpeed 5.X RT service telephone line identified and installed for InSite? (Electrical, mechanical, etc.)

Table 1-6 Miscellaneous Tasks

GE	Cust	Other
Y N	Y N	<i>Must be completed before the scheduled delivery date</i>
	☐ ☐	Arrangements made in the schedule to allow adequate time for remodeling, if required (such as wall, floor, or ceiling repair work, painting, other cosmetic finishes)
	☐ ☐	Have arrangements been made to clean the floor <i>after</i> equipment removal and <i>prior</i> to the installation of the new equipment?
☐ ☐	☐ ☐	Is de-installation of existing equipment required? No ____ Yes ____ Removal date _____

Table 1-6 Miscellaneous Tasks (Continued)

GE	Cust	Other
Y N	Y N	<i>Must be completed before the scheduled delivery date</i>
☐ ☐	☐ ☐	Is there a trade-in of existing equipment? No __ Yes __ GoldSeal _____
☐ ☐	☐ ☐	Delivery route identified and verified with the proper hospital personnel? No__ Yes __ Elevators and doors checked for size and weight constraints? No__ Yes __
☐ ☐	☐ ☐	Have appropriate arrangements been made with traffic for delivery?No__ Yes __
☐ ☐	☐ ☐	Will acceptance/performance testing or bio-medical testing be required?No__ Yes __ Date: _____
☐ ☐	☐ ☐	Are trash and/or recycling bins available for the removal of papers, boxes, etc. during the installation? No__ Yes __

Section 4.0

Medical Electrical Equipment for EMC

4.1 General Scope

This LightSpeed RT System complies with IEC60601-1-2 Edition 2 EMC standard for medical electrical equipment.

The LightSpeed RT System is suitable to be used in the electromagnetic environment, as per the limits & recommendations described in the tables hereafter:

- Emission Compliance level & limits (see [Table 1-7](#))
- Immunity Compliance level & recommendations to maintain equipment clinical utility (see [Table 1-8](#), [Table 1-9](#), and [Table 1-10](#)).

Note: This system complies with above mentioned EMC standard when used with supplied cables up to maximum lengths referenced in the MIS MAPS or system cable interconnect diagrams.

4.2 Electromagnetic Emission

Table 1-7 Emission Declaration

EMC Emissions Guidance & Declaration for LightSpeed RT System		
The LightSpeed RT system is intended for use in the electromagnetic environment specified below. The customer or the user of the LightSpeed RT system should assure that it is used in such an environment.		
Emissions Test	Compliance	Electromagnetic Environment Guidance
RF emissions CISPR 11	Group 1	The LightSpeed RT system uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class A	
Harmonic emissions IEC 61000-3-2	Not applicable	The LightSpeed RT system is suitable for use in all establishments other than domestic and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations/ flicker emissions IEC 61000-3-2	Not applicable	

4.3 Electromagnetic Immunity

Table 1-8 Immunity Declaration

EMC Immunity Guidance & Declaration for LightSpeed RT System			
The LightSpeed RT system is intended for use in the electromagnetic environment specified below. The customer or the user of the LightSpeed RT system should assure that it is used in such an environment.			
Immunity Test	IEC 60601-1-2 Test Level	Compliance Level	Electromagnetic Environment Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete, or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines ± 1 kV for input/output lines	± 2 kV for power supply lines ± 1 kV for input/output lines	Mains power quality should be that of typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV line-line ± 2 kV line-earth	± 1 kV line-line ± 2 kV line-earth	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	< 5 % U_T (> 95% dip in U_T) for 5 sec	< 5 % U_T (> 95% dip in U_T) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the LightSpeed RT system requires continued operation during power mains interruptions, it is recommended that the LightSpeed RT system be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note : U_T is the a.c. mains voltage prior to application of the test level.			

Table 1-9 Immunity Declaration con't

EMC Immunity Guidance & Declaration for LightSpeed RT System			
The LightSpeed RT system is intended for use in the electromagnetic environment specified below. The customer or the user of the LightSpeed RT system should assure that it is used in such an environment.			
Immunity Test	IEC 60601-1-2 Test Level	Compliance Level	Electromagnetic Environment Guidance
Radiated RF IEC 61000-4-3 (alternative method: IEC 61000-4-21)	3 V/m 80 MHz to 2.5 GHz	3 V/m 80 MHz to 2.5 GHz	$d = \left[\frac{3.5}{3} \right] \sqrt{P}$ 80 MHz to 800 MHz (see Table 1-10) $d = \left[\frac{7}{3} \right] \sqrt{P}$ 800 MHz to 2.5GHz (see Table 1-10) where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b. Interference may occur in the vicinity of equipment marked with the following symbol: 
a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast, and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the LightSpeed RT System is used exceeds the applicable RF compliance level above, the LightSpeed RT System should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the LightSpeed RT System. b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m. Note : These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.			

Table 1-9 Immunity Declaration con't

EMC Immunity Guidance & Declaration for LightSpeed RT System			
The LightSpeed RT system is intended for use in the electromagnetic environment specified below. The customer or the user of the LightSpeed RT system should assure that it is used in such an environment.			
Conducted RF IEC 61000-4-6	3 V_{RMS} 150kHz to 80 MHz	3 V 150kHz to 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of the LightSpeed RT System, including cables, than the recommended separation distance calculated from the equation appropriate for the frequency of the transmitter. Recommended Separation Distance $d = \left[\frac{3.5}{3} \right] \sqrt{P}$ (see Table 1-10)
Radiated RF IEC 61000-4-3 (alternative method: IEC 61000-4-21)	3 V/m 80 MHz to 2.5 GHz	3 V/m 80 MHz to 2.5 GHz	$d = \left[\frac{3.5}{3} \right] \sqrt{P}$ 80 MHz to 800 MHz (see Table 1-10) $d = \left[\frac{7}{3} \right] \sqrt{P}$ 800 MHz to 2.5GHz (see Table 1-10) where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b. Interference may occur in the vicinity of equipment marked with the following symbol: 
a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast, and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the LightSpeed RT System is used exceeds the applicable RF compliance level above, the LightSpeed RT System should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the LightSpeed RT System. b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m. Note : These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.			

Table 1-10 Separation Distances

Recommended separation distances between portable and mobile RF communications equipment and the LightSpeed RT System.			
The LightSpeed RT System is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the LightSpeed RT System can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the LightSpeed RT System as recommended below, according to the maximum output power of the communications equipment.			
Rated Maximum Output Power (P) of Transmitter Watts (W)	Separation distance according to frequency of transmitter		
	150 kHz to 80 MHz	80 MHz to 800MHz	800 MHz to 2.5 GHz
	$d = \left[\frac{3.5}{3} \right] \sqrt{P}$	$d = \left[\frac{3.5}{3} \right] \sqrt{P}$	$d = \left[\frac{7}{3} \right] \sqrt{P}$
	Separation Distance meters	Separation Distance meters	Separation Distance meters
0.01	0.12	0.12	0.23
0.1	0.37	0.37	0.74
1	1.17	1.17	2.33
10	3.69	3.69	7.38
100	11.7	11.7	23.3
For transmitters rated at a maximum output power not listed above, the separation distance can be estimated using the equation in the corresponding column, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer. Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies. Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.			

4.3.1 Limitations Management :

Adhering to the distance separation recommended in [Table 1-10](#), between 150kHz & 2.5GHz, will reduce disturbances recorded at the image level but may not eliminate all disturbances. However, when installed and operated as specified herein, the system will maintain its essential performance by continuing to acquire, display, and store diagnostic quality images safely.

(*) For example, a 1W mobile phone (800MHz to 2.5GHz carrier frequency) shall be put 2.3 meters apart from the LightSpeed RT System (in order to avoid image interference risks).

4.4 Use Limitation :

4.4.1 External components

The use of accessories, transducers, and cables other than those specified may result in degraded ELECTROMATHNETIC COMPATIBILITY of the LightSpeed RT System.

4.5 Installation Requirements & Environment Control :

In order to minimize interference risks, the following requirements shall apply.

4.5.1 Cable Shielding & Grounding

All interconnect cables to peripheral devices must be shielded and properly grounded. Use of cables not properly shielded and grounded may result in the equipment causing radio frequency interference.

4.5.2 This product complies the radiated emission as per CISPR11 Group 1 Class A standard limits

The LightSpeed RT System is predominantly intended for use, in non-domestic environments, and not directly connected to the Public Mains Network. The LightSpeed RT System is predominantly intended for use (e.g. in hospitals) with a dedicated supply system, and with a X-ray shielded room. In case of using in a domestic environment (e.g. doctor's offices), in order to avoid interferences, it is recommended to use a separated AC power distribution panel & line with a X-ray shielded room.

4.5.3 Subsystem & Accessories Power Supply Distribution

All components, accessories subsystems, systems which are electrically connected to the LightSpeed RT System, must have all AC power supplied by the same power distribution panel & line.

4.5.4 Stacked Components & Equipment

The LightSpeed RT System should not be used adjacent to or stacked with other equipment; if adjacent or stacked use is necessary, the LightSpeed RT System should be observed in order to verify normal operation in the configuration in which it will be used.

4.5.5 Low Frequency Magnetic Field

In case of a digital LightSpeed RT System, the Gantry (digital detector) shall be apart 1 meter from the generator cabinet, and 1 meter apart from the analog (CRT) monitors. These distance specifications will minimize the low frequency magnetic field interference risk.

4.5.6 Static Magnetic Field Limits

In order to avoid interference on the LightSpeed RT System system, static field limits from the surrounding environment are specified.

Static field is specified less than <1 Gauss in Examination room, and in the Control Area.

Static field is specified less than <3 Gauss in the Technical Room.

4.5.7 Electrostatic Discharge Environment & Recommendations

In order to reduce electrostatic discharge interference, install a charge dissipative floor material to avoid electrostatic charge buildup.

The relative humidity shall be at least 30 percent.

The dissipative material shall be connected to the system ground reference, if applicable.

Chapter 2

Pre-Installation Overview



- NOTICE** Prior to beginning the pre-installation process for a LightSpeed RT system, it is highly recommended that the user review the Pre-Installation Video (DVD).
- Before a LightSpeed RT system can be installed, all pre-installation requirements must be complete.
- Chapter 4, Sections [8.0](#) and [9.0](#) Structural Requirements
 - Chapter 4, Section [10.0](#) Radiation Protection
 - Broadband standard
 - Site Ready Visit
 - [Chapter 5, Environmental Conditions](#), Sections [1.0](#) & [2.0](#) HVAC Requirements
 - [Chapter 8, Power Requirements](#) (Site Power Audit Required)
 - [Chapter 9, Interconnection Data](#)

Site-specific items must be verified before the installation can begin.

Section 1.0

Dust/Dirt Contamination

The LightSpeed RT system (consisting of: Console, PDU, Table and Gantry) is highly susceptible to airborne contaminants, especially concrete and drywall dust. Due to the possibility of contamination, these systems should NEVER be installed in a construction site.



- NOTICE** Any site with unfinished floors, walls or ceilings is considered a construction site, and is not suitable for system installation.

Section 2.0

Chemical Contamination

Wet film processors must never be installed in the same room as the scanner, due to the possibility of chemical contamination of LightSpeed Series components. Such chemicals can contribute to increased equipment failures, increased system downtime, and decreased reliability. Film processor equipment installation must meet the manufacturer's requirements (e.g. ventilation specifications) and all applicable national and local codes. Also, consideration's should be given to the location of this equipment and chemical fumes relative to human contact as it relates to locating this equipment and chemicals in the control room.

Section 3.0

Walls, Ceiling, and Floor

All walls, ceiling, and flooring must be completed before installation can begin.
LightSpeed RT scanners can only be installed on a 4" concrete floor.

Section 4.0 Broad-band

For information on Broad-band requirements, refer to Chapter 4, [Section 9.0, on page 74](#).

Section 5.0 Phone Line (for optional modem)

Two phone lines must be installed at or near the console and be operational prior to installation.

- 1-Analog line (for modem use)
- 1-Voice line

Section 6.0 Review

LightSpeed RT systems use adjustable leveling pads to support the gantry and table. The gantry has four (4) primary leveling pads. The table has five (5) pads used for leveling it.

Using the GE print to establish the room layout, make sure all the operating and service clearances shown on the print are observed. Using the template (P/N 2389013) shipped with the system, locate the anchor holes. Make sure they clear structural interferences in the floor.

Clean the area. Free the mounting surface of any material that may interfere with the positioning and leveling of the system.

- 1.) Lay out the 2 floor templates.
- 2.) Start with the Gantry template—align per the GE print.
- 3.) Place the Gantry template over the top of the Table template. Align the scan and table center-lines and secure the templates to the floor. Make sure there are no potential clearance issues.
- 4.) Check the level of the floor (See [Figure 2-1](#)) across the templates.

Note: Tiles (or other resilient flooring) around all holes will be cut during the installation process.

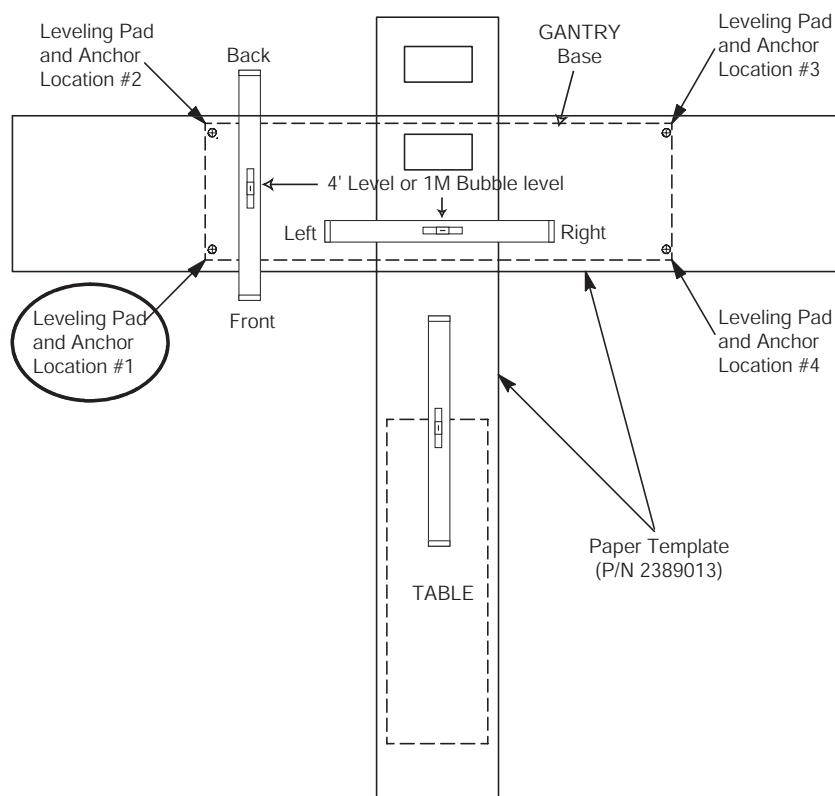


Figure 2-1 Hole Locations

FLOOR LEVELNESS SPECIFICATION

- 3.16 mm (1/8 in.) over 3048 mm (10 ft.)

This should be measured on the template over the table/gantry area, as shown in [Figure 2-1](#), above..

This Page Intentionally Blank.

Chapter 3

System Catalog

Section 1.0

Option Catalog Numbers

The following is a list of system options requiring site planning work for the LightSpeed RT system. Contact your local GE Healthcare Sales representative for a complete list of all system options or visit us at <http://www.gehealthcare.com>. Refer to the instruction manuals supplied with specific options for respective details.

CATALOG NUMBER	OPTION DESCRIPTION
International Dolly Set B7850LD	For International customers, if dollies are required.
Rear Cable Cover B7850RC	Optional rear cable cover for LightSpeed RT upgrades
Gammex Lasers E8505LF	CT Simulation Laser Systems
LAP Lasers E8505P	CT Simulation Laser Systems
Exact Couch E8505MA	Radiation Therapy Planning (RTP) Exact Couch for Non-NP Systems. Flat panel table top provides a flat surface for accurate and reproducible patient positioning.
Remote Color Monitors B7710WM B7530RC	LightSpeed RT LCD monitor. Remote 20 inch diagonal color monitor.
Performance Network Kits K9000L	6 Node, 10/100 Mbit Auto Sensing
Switched Network Kit B7500PM	ConnectPro Option provides a direct interface to HIS/RIS
Bar Code Reader B7540RB	LightSpeed RT Bar Code Reader
Seismic Kit R4390JC	Earthquake zones
Long Cable Set B7816PH	

Table 3-1 System Options Catalog (Part) Numbers

Section 2.0

Base Scanner System

2.1 Application

The CT scanner system includes hardware and software to support patient data acquisition and image analysis for whole-body computed tomography.

2.2 Configuration

The base scanner system is configured as shown. All scan and analysis functions are controlled from the operator's console (not shown).

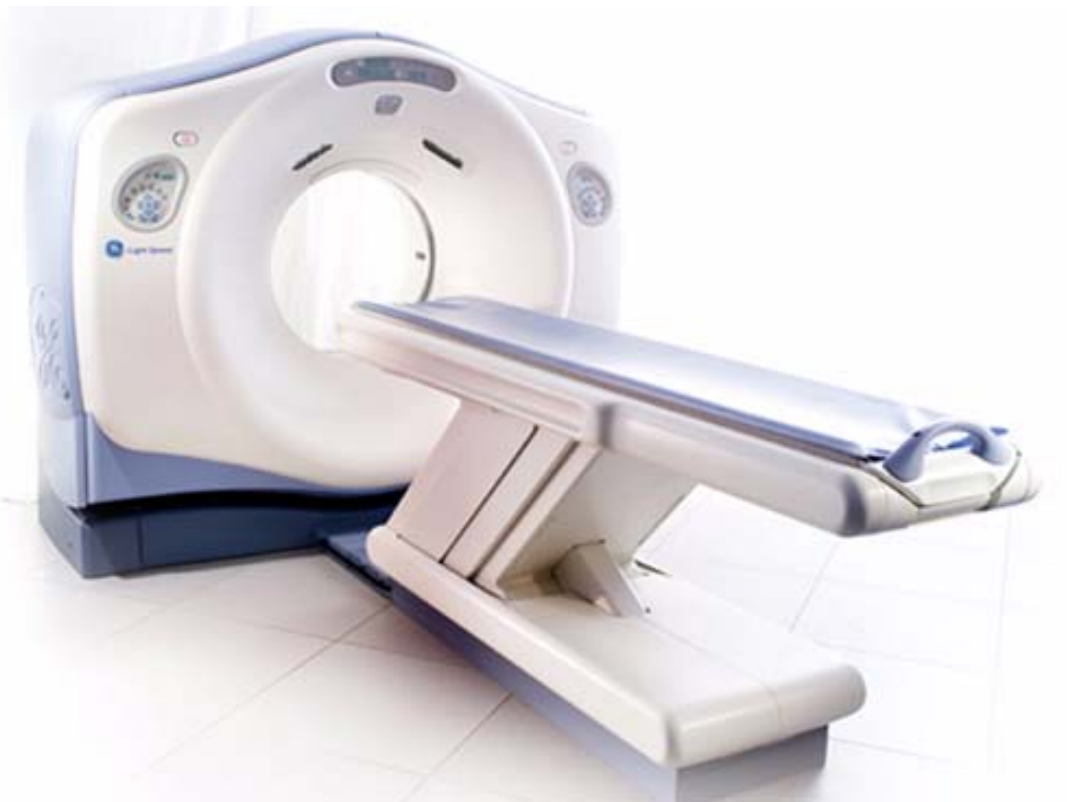


Figure 3-1 Base Scanner System

Chapter 4

Room Planning

Section 1.0 Required Systems Clearances



NOTICE

If your system was installed before December 2006, Limited Access is not available for your system. You must use the Recommended Room Size to reinstall it.

Consult your local GE Sales and Service Representative about your specific needs.

System Configuration	Recommended Scan Room Size	Typical Scan Room Size	Recommended Ctrl Room Size	Clearance Requirements
HP Table	4420 x 6096 mm (14 ft. 6 in. x 20 ft.)	4267 x 5867 mm (14 ft. x 19 ft. 3 in.)	2743 x 4420 mm (9 ft. x 14 ft. 6 in.)	see Figure 4-1

Limited Access Rooms (width only):

Minimum Width: 3708 mm (12 ft. 2 in.) 356 mm (14 in.) cover to wall

Average Width: 4064 mm (13 ft. 4 in.) 711 mm (28 in.) cover to wall

Additional component dimensions are available in [Figure 4-8](#) through [Figure 4-12](#) of this document. Consult your local General Electric Project Manager of Installation for your appropriate room specifications. For equipment clearance requirements, refer to [Section 2.0](#). Remember, sufficient Regulatory and Service clearances must be maintained around equipment for full operation, service and safety. Cable length is an important consideration in room layout. The LightSpeed RT system is shipped with standard short length cables, with a set of longer cables (B7816PH) available as an option. See the electrical page of the GE print for your specific requirements.

Note, also, that where possible, the cables should enter the Gantry from the rear, utilizing the rear cable cover assembly. Alternate cable entry is possible at the center of the Gantry (refer to the Installation template).

- Excess cable length cannot be stored behind the Operator Console or PDU.
- Long cable must not be cut or shortened.
- Excess cable may be stored in conduits, a cable storage box if present, or the floor duct, provided sufficient space is available. Observe the fill rate for each option. If there are questions regarding local electrical or building codes, consult the project electrical contractor or electrician.
- All NEC 70-E Electrical Regulations regarding conduit or duct fill must be observed.

Section 2.0 Regulatory and Service Clearances

2.1 Regulatory Clearances

MINIMUM CLEARANCES UNDER U.S. FEDERAL REGULATIONS AND NATIONAL STANDARDS:

29 CFR 1910 (OSHA),

NFPA 70E (STANDARD FOR ELECTRICAL SAFETY IN THE WORKPLACE)

NFPA 101 (LIFE SAFETY CODE):

Figure 4-1 is a map of clearance requirements for U.S. regulatory compliance. See clearance tables on the following pages for detailed dimensional clearances. Please note all systems installed in the United States must comply with all Federal and local regulations. For installations outside the United States, country-specific or other local regulatory clearance requirements must be met.

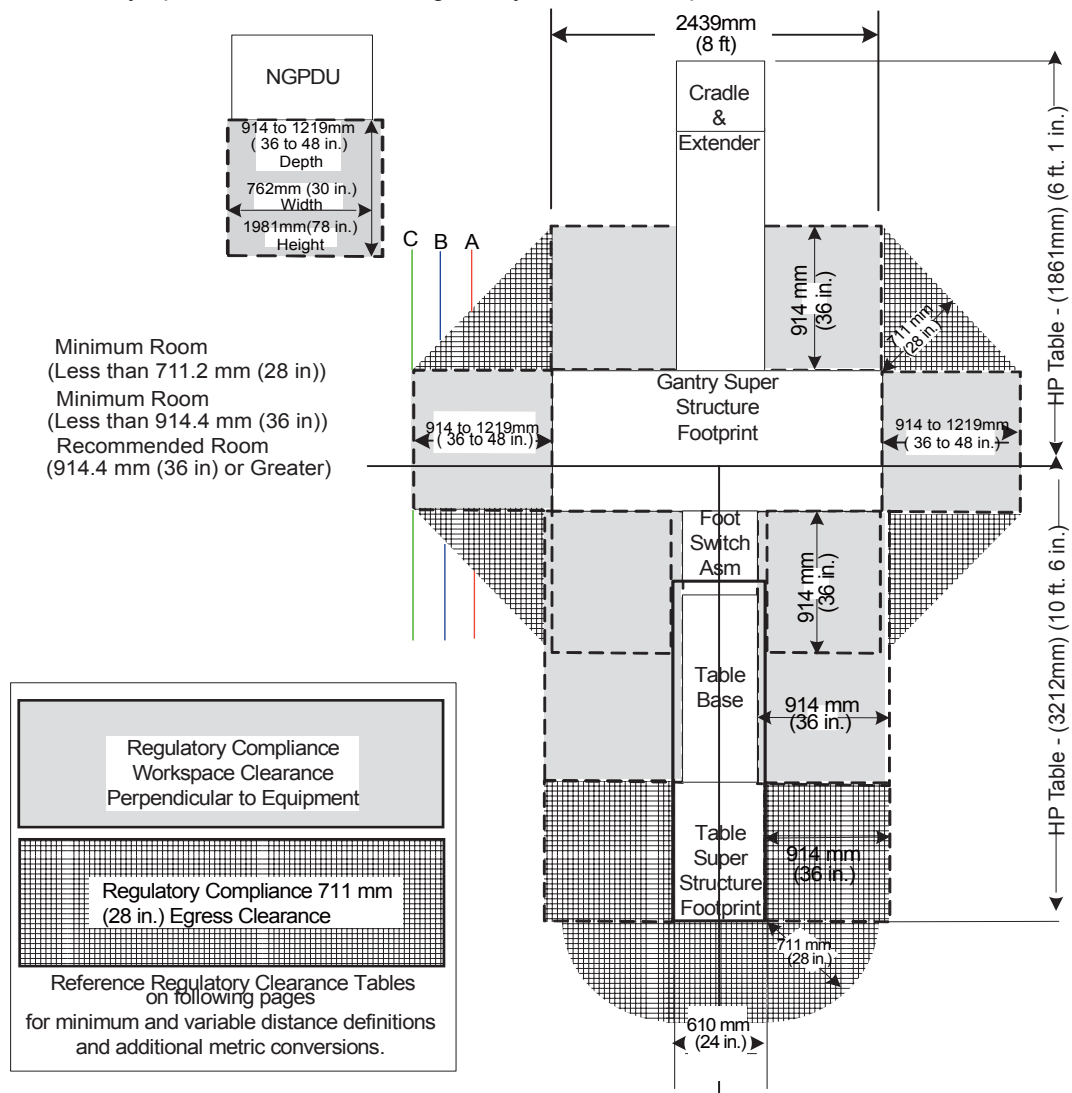


Figure 4-1 Regulatory Clearance Requirements for LightSpeed RT System Configuration

Note: See [Section 4.0](#) for Service Clearances

2.1.1 Regulated Minimum Working Clearance by Major Subsystem

- Requirements apply to equipment operating at 600V or less, where examination, adjustment, servicing, or maintenance is likely to be performed while live parts are exposed.
- Direction of Service Access is defined as perpendicular to the surface of the equipment being serviced.
- Required regulatory clearance distances must be maintained and may not be used for storage. This includes normal system operation as well as service inspection or maintenance..

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access: front of console	914 mm (36 in.)	There are no exposed live part hazards with the cover in place. If the console is placed under a counter, the front edge of the console must be even with the vertical edge of the console workspace. Note: This component is typically serviced from the front with access to the rear.
Service access width: Front of console	762 mm (30 in.)	This is the width of the workspace in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.
Head clearance	1981.2 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981.2 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 4-1 Console Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Front of NGPDU)	914.4 mm (36 in.)*	There are no exposed live part hazards with the cover in place. This component is typically serviced from the front with access to the rear. *If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) is required on both sides of the workspace with the operator between. *If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.
Service Access Width (Left-Right of workspace)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.

Table 4-2 NGPDU Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Head Clearance	1981 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 4-2 NGPDU Subsystem

- For the gantry and table, distances are measured from the enclosure, not the finish covers.

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (All Sides)	914 mm (36 in.)	If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) on both sides of workspace with the operator between is required. If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.
Service Access Width (Left-Right of workspace)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.

Table 4-3 Gantry Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Table Head or Foot)	914 mm (36 in.)	There are no exposed live parts hazards with the cover in place. This component is typically serviced from all four sides. this is the width of the workspace on each side of the equipment. A minimum of 914.4 mm (36 in.), or the width of the equipment, whichever is greater, is required.
Direction of Service Access (Table Sides)	914 mm (36 in.)*	*This distance can be reduced to 711 mm (28 in.) provided a written and signed approval is obtained by the local team from the local AHJ (Authority Having Jurisdiction). The signed document must be on file with GE.

Table 4-4 Table Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Table Foot)	711 mm (28 in.)	For the front gantry cover removal, a minimum of 457 mm (18 in.) is allowed only if an unobstructed egress space of 711 mm (28 in.) is maintained around the equipment for room exit. This also means no trip hazards exist along the path of egress.
Service Access Width (Left-Right of workspace)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.

Table 4-4 Table Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Front of UPS)	914.4 mm (36 in.)*	There are no exposed live part hazards with the cover in place. This component is typically serviced from the front with access to the rear. *If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) is required on both sides of the workspace with the operator between. *If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.
Service Access Width (Right side and length of UPS)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.
Head Clearance	1981 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 4-5 UPS Subsystem

Work Space Requirement	Minimum Clear Space	Additional Conditions
Direction of Service Access (Front of A1 Disconnect)	914.4 mm (36 in.)*	There are no exposed live part hazards with the cover in place. This component is typically serviced from the front with access to the rear. *If exposed live parts of 151 - 600 volts are present, 1219 mm (48 in.) is required on both sides of the workspace with the operator between. *If the opposite wall is grounded and exposed live parts of 151 - 600 volts are present, 1067 mm (42 in.) is required.
Service Access Width (Right side and length of A1 Disconnect)	762 mm (30 in.)	This is the width of the working space in front of the equipment. A minimum of 762 mm (30 in.) or the width of the equipment, whichever is greater, is required.
Head Clearance	1981 mm (78 in.)	This is the height of the workspace measured from the floor at the front edge of the equipment to the ceiling or overhead obstruction(s). A minimum of 1981 mm (78 in.) or the height of the equipment, whichever is greater, is required.

Table 4-6 A1 Disconnect Subsystem

2.1.2 Terms and Definitions

EGRESS

The path of exit from within any room. U.S. regulatory requires a minimum of 711.2 mm (28 in.) of continuous and unobstructed space including trip hazards along the path of exit.

WORK SPACE

This is the dimensional box required for safe inspection or service of energized equipment. It consists of depth, width, and height. The depth dimension is measured perpendicular to the direction of access. U.S regulation is minimum of 914 mm (36 in.). Additional conditions can increase the minimum requirement. FCT defines this as the envelope of the component superstructure. For the NGPDU it is with the front panel removed. For the gantry and table, it is with the patient or external covers removed.

SERVICE ACCESS WIDTH

This is the width of the working space in front of the equipment, a minimum of 762 mm (30 in.), or the width of the equipment whichever is greater.

HEAD CLEARANCE

This is the height dimension of "Work Space". The height of the workspace measured from floor at the front edge of equipment to ceiling or overhead obstruction(s), 1981 mm (78 in.) or height of equipment, which ever is greater.

GROUNDING WALL

Any wall that can be electrically conductive to earth ground. Masonry, concrete, or tile, are considered conductive. Additional commonly found aspects of a wall should also be considered as grounded. This is not an all-inclusive list:

- Medical Gas ports
- Metal door and window frames
- Water sources and metallic sink structures
- Metallic wall mounted cabinets
- A1 disconnect panel
- Equipment Emergency Off panels
- Industrial equipment such as air conditioners and vents
- Expansion joints

The following are not considered as grounded elements of a common wall:

- Standard wall outlet
- Light switches
- Telephones
- Communication wall jacks

MINIMUM

The lowest limit permitted by law or other authority.

DIMENSIONS AND CLEARANCES

Consisting of, or representing the lowest possible amount of degree for freedom permissible for equipment siting. This relationship must meet all safety, service, and regulatory requirements to be acceptable.

PRE-INSTALLATION ESCALATION

Process to consult with CT Engineering, the Design Center or EHS regarding pre-installation issues related to your siting concerns.

Section 3.0

Additional Regulatory Clearance Information

3.1 Minimum Room Size (Limited Access)

The CT Gantry Left Side Limited Access Initiative provides the capability to reduce the minimum room size for CT Systems while still meeting all installation requirements and specifications. This adds left side flexibility, allowing the CT system to be sited in rooms with widths 558.8 mm (22 in.) smaller than the current minimum room width. Left-side access and egress may be restricted. Refer to your site's installation print for your room's detail.

If you are using the square meters (square feet) to determine compliance, remember that the front cover clearance is wider than the regulatory clearance along the table length, and that the cover park position is behind the table in the home position. The cover is removed on tilting dollies and can be moved side-to-side to reach the park position of 457 mm (18 in.), if this is not an egress route.

Wall duct and conduit on walls within the regulatory clearance areas shall be 1067mm (42 in.), measured from the covers to the obstruction. Servicing of the CT System can be safely performed within the regulatory envelopes, however sufficient space must be maintained to remove system covers, and replace large system components. To achieve this clearance for the gantry, clear space must be available to maneuver the gantry covers mounted on the service dollies. Surface floor raceway cannot be used in the egress route areas. OSHA ramps are available. The FE lifting the rear or front cover to avoid floor obstructions is not an EHS-approved service procedure.

One Service Engineer shall be able to accomplish all service component replace tasks listed without the need for special tools or equipment, such as a tube change, detector change, and HV tank.

3.1.1 Regulatory Caution

Site prints are required for all system installations including relocation and moves. CT room layout, as shown on your site print, shall meet all regulatory requirements as described in the installation manual. Additional room components, such as cabinets, reduce room size. Equipment not shown on the site print may void the caution statement, making the room non-compliant. Actual site measurements before installation will be taken to determine room size and compliance.

3.1.2 Egress Clearance

Egress requires a clear, unobstructed route out of the room, either around the back of the gantry or around the back of the table. If your egress route is not around the back of the table, maintain 457 mm (18 in.) of clearance between the back of the table, with a continuous width of 3200 mm (126 in.), 1600 mm (63 in.) on each side of the table center line, on each side to any obstruction so that the front cover can be removed. Refer to the Pre-Installation manual for more details on service clearances.

Exceptions

Rooms smaller than 3708mm x 5867mm (12 ft. 2 in. x 19 ft. 3 in.) require construction to meet the minimum requirements. The design center or your GE PMI may have additional recommendations for your room size.

3.1.3 Operational Caution

In a minimum room layout 355.6 mm - 685.8 mm (14 in. - 27 in.), the customer should consider workflow, customer access for patient care, and critical-care operations space requirements. Additionally, there may be limited equipment access on the gantry left side when loading patients or when positioning patient equipment in the room between the gantry and the wall. Detailed customer installation tasks are detailed in the product Pre-Installation manual, Chapters 1-4.

3.1.4 System Specifications (LightSpeed RT Systems)

System Configuration	Recommended Scan Room Size	Typical Scan Room Size	Minimum Room Size
HP Table	4420 x 6096 mm (14 ft. 6 in. x 20 ft.)	4267 x 5867 mm ¹ (14 ft. x 19 ft. 3 in.)	3708 x 5867 mm ² (12 ft. 2 in. x 19 ft. 3 in.)

¹Same Regulatory requirements apply

²Same Regulatory requirements apply, with the addition of no energized left side service.

3.1.4.1 Recommended Room Size

The "Recommended room size" configuration offers the most flexibility for future upgrades. It has sufficient workspace and space to add millwork, while meeting all regulatory requirements. This room is compatible with most two-step future installations.


3.1.4.2 Typical Room Size

The "Typical room size" configuration allows for some future upgrades. It has sufficient workspace, but limited space to add millwork and meet all regulatory requirements. This room may be compatible with some two-step future installations.

3.1.4.3 Minimum Room Size

The "minimum room size" configuration allows for no future upgrades. It has limited workspace and no in-room millwork, but meets all regulatory requirements. This room is not compatible with two-step future installations.


3.2 How to Measure

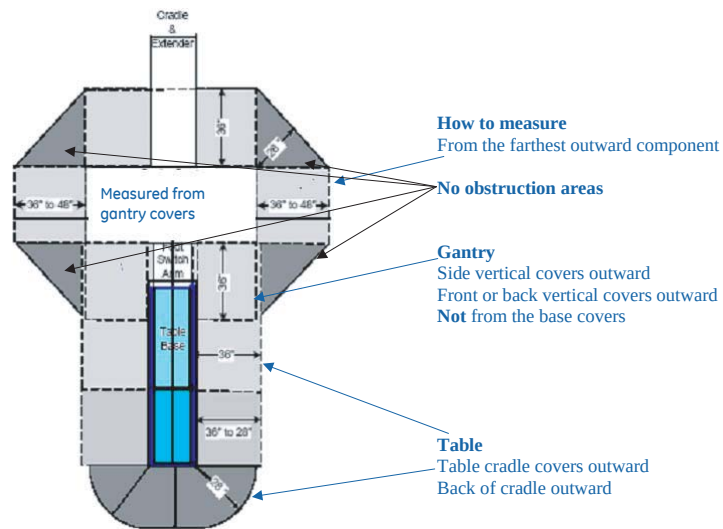

3.3 Minimum Room Size & Requirement Layouts


Room A - Less than 711 mm (28 in.), but greater than 256 mm (14 in.), measured from the covers to the left sidewall.

In this configuration service, egress and workspace are compromised around the gantry's left side.

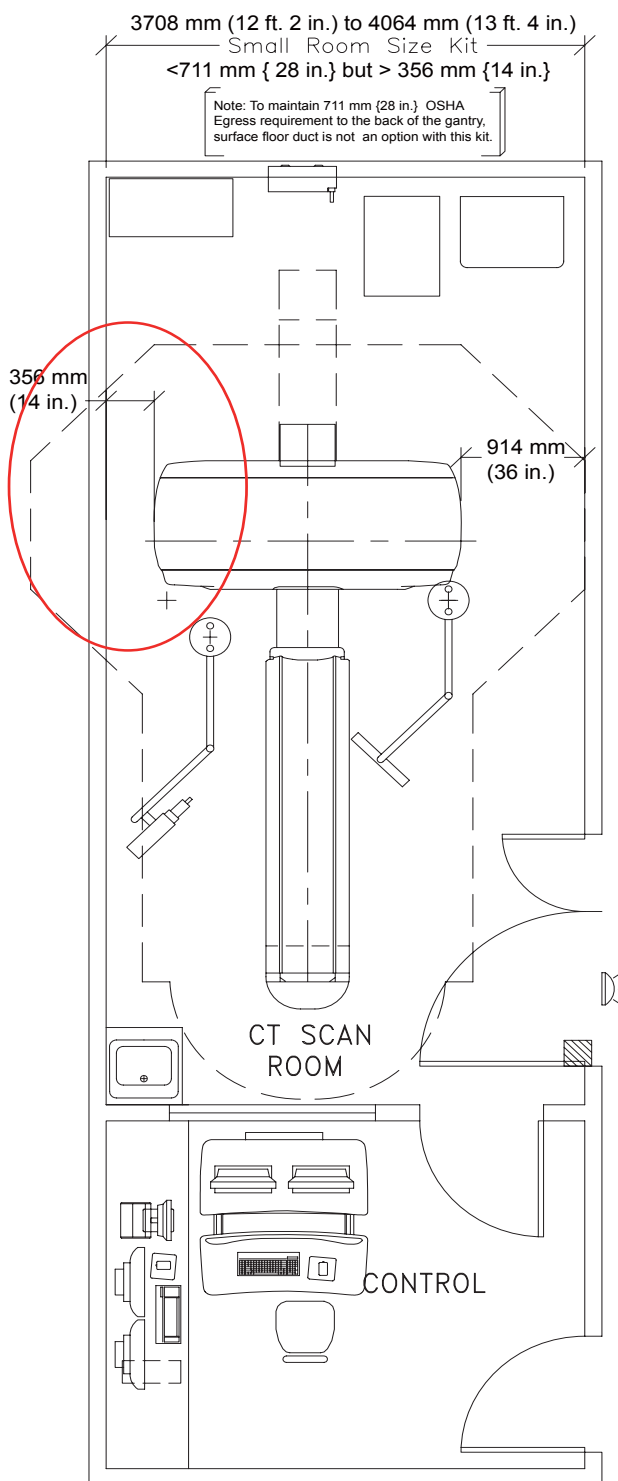

Room B - Less than 914 mm (36 in.) but greater than 711 mm (28 in.) measured from the covers to the left sidewall.

In this configuration service, egress, and workspace are acceptable around the gantry.

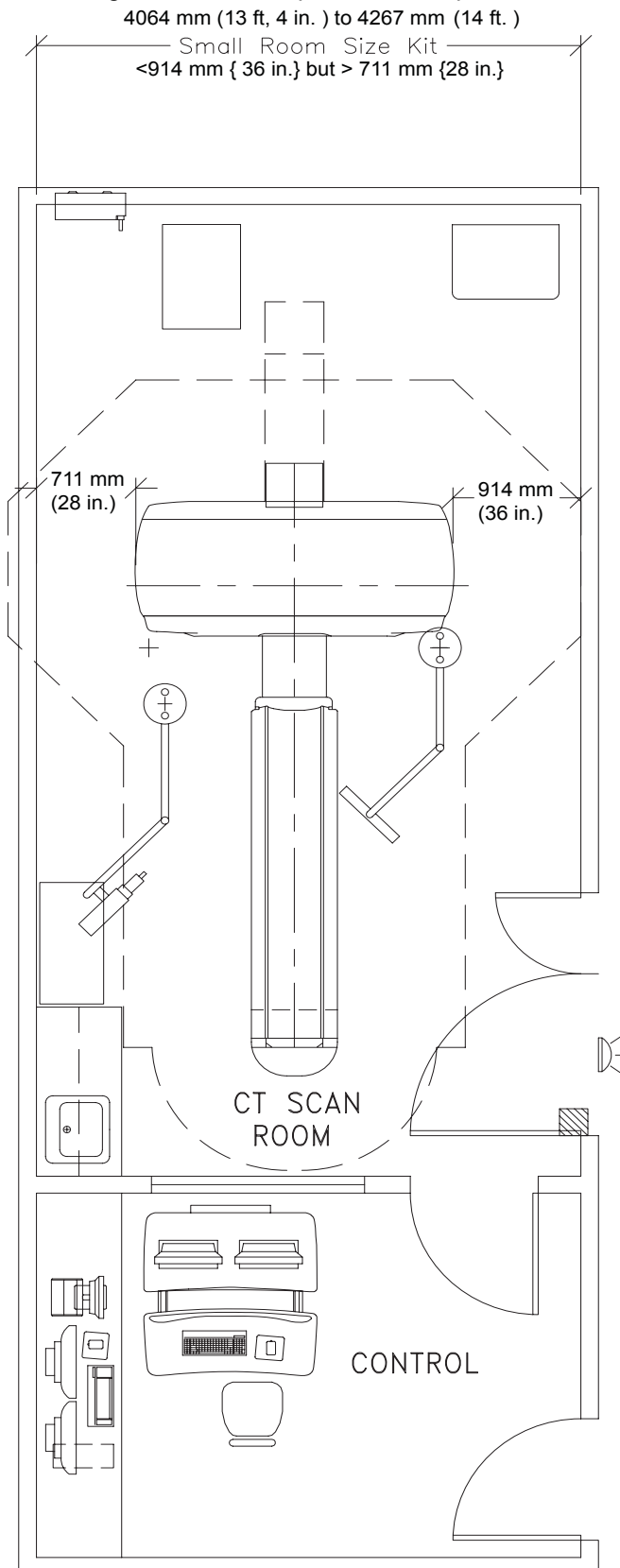

3.4 Recommended and Typical Room Size & Requirement Layouts

Figure 4-2 and Figure 4-3 show Recommended and typical room layouts, with and without in-room cabinets. You need to be aware of locations for me gas, surface ductwork, or other items that make a grounded wall.

Note: Your room layout may meet the Recommended or typical room requirements but look different than Figure 4-2 or Figure 4-3. Contact your sales person to have a detail room layout completed for your site.

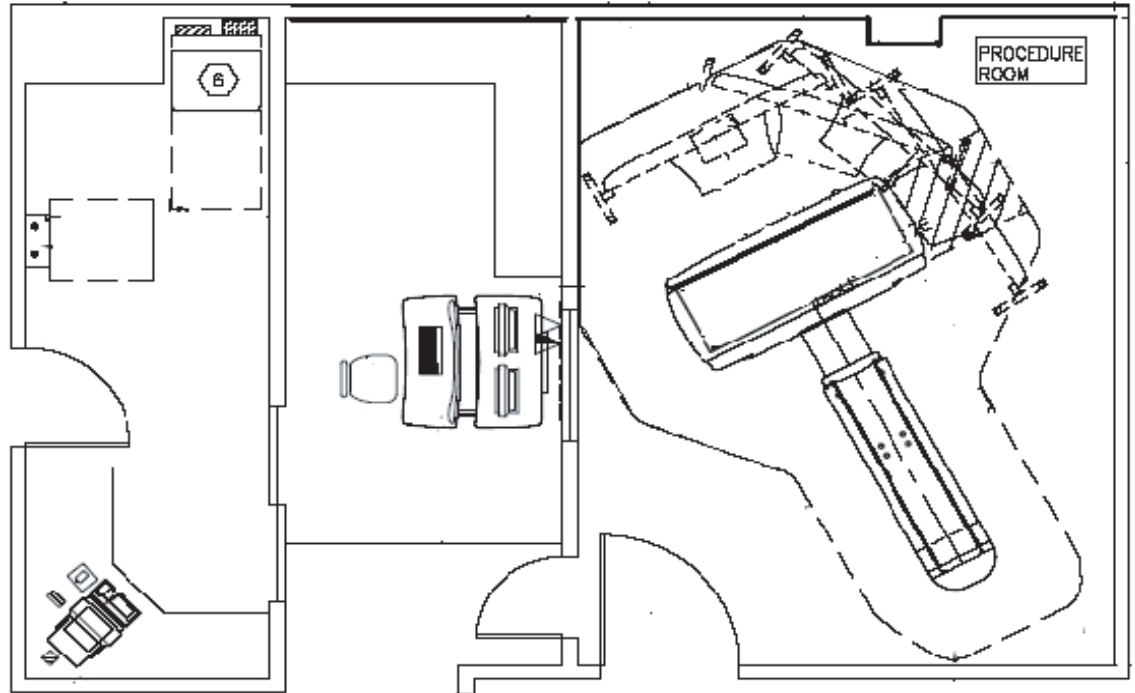


Figure 4-2 Room Layout with cabinets

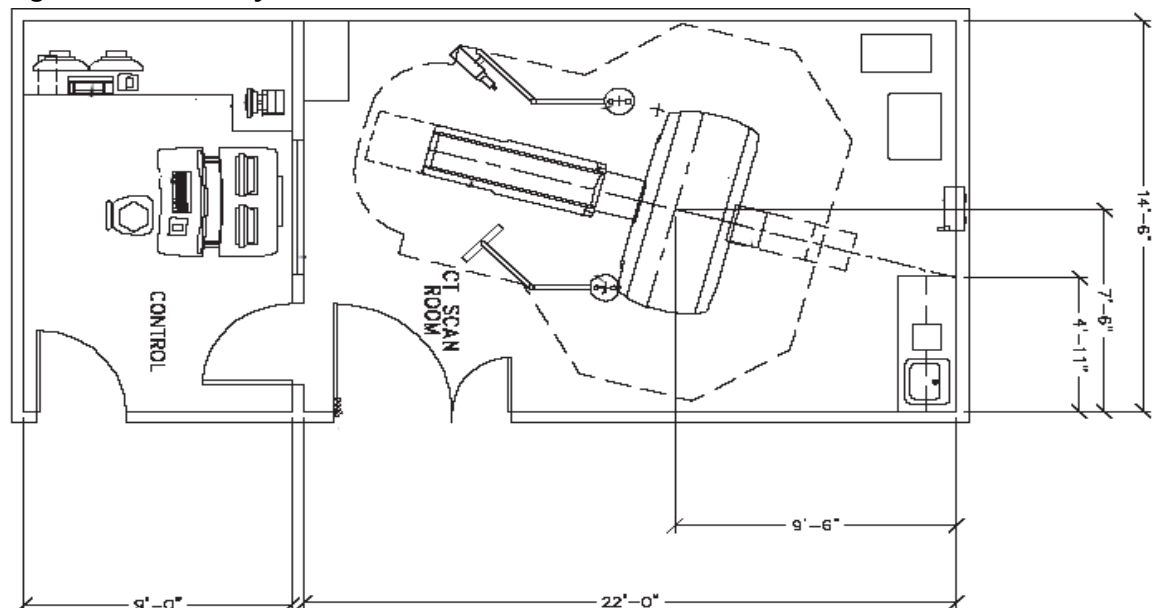


Figure 4-3 Room Layout without cabinets

Section 4.0 Service Clearance

Servicing of the CT System can be safely performed within the regulatory envelopes defined in [Section 2.1](#), however sufficient space must be maintained to remove the covers from the system. To achieve this clearance for the gantry, clear space must be available to maneuver the gantry covers mounted on the service dollies. One Service Engineer can accomplish this. If such space is not available then two Service Engineers will be required for gantry cover handling. This must be clearly communicated with the Customer.

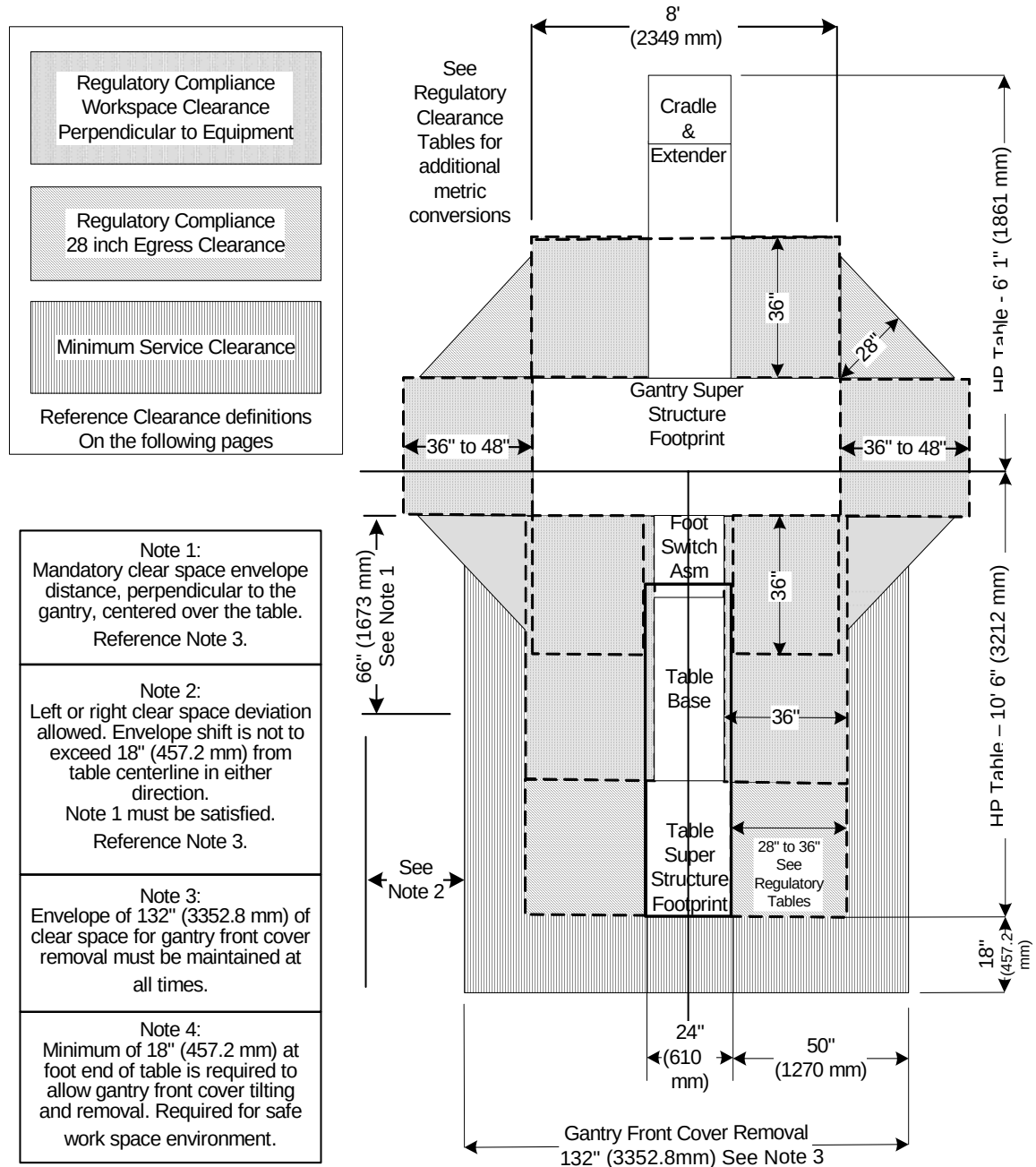


Figure 4-4 Minimum Service Clearance

4.1 Service Clearances for Single Service Engineer

- Gantry front cover removal requires the use of the "Tilting Cover Dollies". These dollies allow the Service Engineer to separate the cover from the Gantry, tilt the cover 90 degrees, roll the cover to the foot end of the Table, and then tilt the cover an additional 90 degrees such that the front cover is now upside down relative to the normal system mounted condition. [Figure 4-4](#) illustrates the minimum clear space required to achieve this operation of 3,200.4 mm (126 in.). The Gantry front cover must be removed to a position that satisfies the minimum regulatory clearances.
- The Gantry rear cover, with service dollies installed, requires a width 2660 mm (104 ft, 7 in.) and a depth of 596 mm (23 ft, 5 in.) of clearance for removal. Sufficient space must be calculated to move the cover either straight back or to a side of the Table to satisfy the minimum service clearance shown in [Figure 4-4](#). This means the rear cover cannot violate the workspace on the rear or either side of the Gantry.
- If Gantry service requires both the front and rear covers be removed, then these covers must be positioned within the room in such a manner as to not violate the egress clearances on any side of the Gantry. This may necessitate removing the covers from within the suite. This should be discussed with the customer and provisions made to accommodate this potential event.
- A single Service Engineer can safely perform servicing of the Table. Sufficient clear space must be available to maintain egress clearances when the Table covers or cradle are removed.
- In your room layout design, service shall have clear and unobstructed access to the gantry tube change area for all major component replacements. These components must be able to reach the service area without lifting or rigging by one service engineer. Major components include:
 - CT X-ray tube in crate
 - High voltage tank(s) in crate
 - Slipping in crate
 - Detector assembly

Be aware of cabinet placement, and how surface floor ducts are used in room configurations.

4.2 Power Distribution Unit (NGPDU) Service Clearance

Positioning of this component must be considered for regulatory compliance as defined in [Section 2.1](#), Regulated Minimum Working Clearance by Major Subsystem.

See Regulatory Tables.

4.3 Console Service Clearance

The console does not present an exposed live parts hazard. However, it is recommended that a minimum working space depth of 1219.2 mm (48 in.) and full width of the console be maintained at all times for service activity. Additionally sufficient space needs to be provided for repositioning of the console and side clearance for rear service access. Egress as well as other service requirement shall be considered when siting the console. See [Figure 4-5](#) for a typical control room layout.

Section 5.0 Recommended Layouts

5.1 Storage Cabinet

A storage cabinet is provided by GE Healthcare to store all supplied service equipment (see [Table 4-7](#) for equipment list). This storage cabinet (457mm D x 914mm W x 1067mm H) (18" D x 36" W x 42" H; ~90 lbs) should be located in the scan room suite area, for easy service access.

ITEM	SIZE	WEIGHT (TOTAL)	
QA Phantom (water filled)	20 cm x 15 cm (7.9 in. x 5.9 in.)	5.5 kg	12 lb
35CM Phantom	35 cm x 7 cm (13.8 in. x 2.8 in.)	8.2 kg	18 lb
48CM Phantom	48 cm x 7 cm (18.9 in. x 2.8 in.)	11.4 kg	25lb
Phantom Holder	25 cm x 25 cm (9.8 in. x 9.8 in.)	3.6 kg	8 lb
FE Box (Purple)	30 cm x 38 cm x 30 cm (11.8 in. x 15 in. x 11.8 in.)	6.8 kg	15 lb
Rear Cover Dollies	158 cm x 82 cm (62.2 in. x 32.3 in.)	11.4 kg	25 lb
Front Cover Dollies	85 cm x 20 cm and 85 cm x 15 cm (33.5 in. x 7.9 in. and 33.5 in. x 5.9 in.)	15.9 kg	35 lb
Install Support Kit (box)	30 cm x 30 cm x 38 cm (11.8 in. x 11.8 in. x 15 in.)	9.1 kg	20 lb
Tube Hoist Kit	77 cm x 8 cm and 38 cm x 15 cm (30.3 in. x 3.1 in. and 15 in. x 5.9 in.)	9.1 kg	20 lb
Balance Weight Kit		33 kg	73 lb

Table 4-7 Equipment to be stored in storage cabinet

5.2 Advantage Workstation (AW)

Refer to Pre-Installation Manual 2111833 and Installation/Service Manual 2111831. For document access, please refer to GE Healthcare's Common Documentation Library web site.

5.3 Control Room Considerations

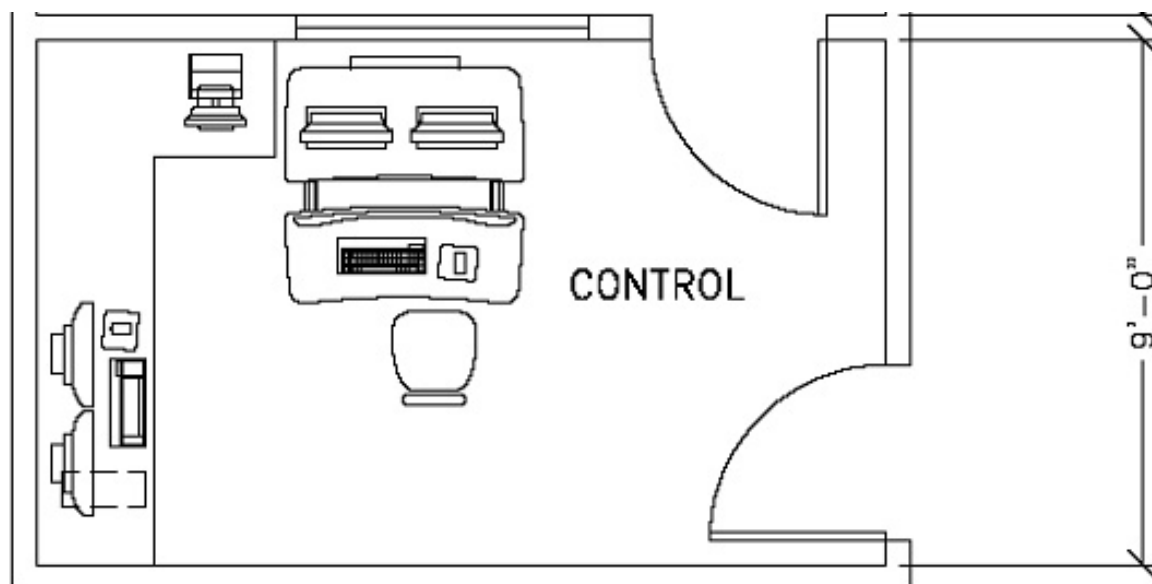


Figure 4-5 Typical Control Room Layout

- The control room must provide a suitable operating environment for the console electronic, and operator working comfort.
- The console cannot be dismantled or have components removed or rearranged in configurations other than as shipped.
- If operationally possible, the monitor desktop and user desktop components may be removed and placed on a counter-top, providing the cable lengths shipped are not altered or changed. The console cabinet then can be remote mounted, provided the cooling and venting requirements are met
- A suitable work area, which is within reach of the Operator's Console, should be provided for placement of the injector control. Injector controls differ in dimensions depending on the brand selected.
- A PACS, workstation, image printer, or filming device are often placed in the console control room area, and sometimes may be directly linked to the console.
- Additional components although linked via network or ethernet cable, are not powered from the CT console.
- Additional room power and network connection should be considered when reviewing the console work space.

5.4 CT Simulation Laser Alignment Systems

CT simulation laser alignment systems are used to reproduce the patient position for simulation and subsequent treatment on a linear accelerator. These lasers can be wall mounted or can be free standing. Work with your installation planner to determine the correct installation solution for your site.

Laser Options: www.gammex.com and www.lap-laser.com

Custom freestanding mounting structures are available for rooms that have no suitable wall or ceiling mounting space. Stanchion systems can provide flexibility in room design for placing alignment devices in rooms that may not be suitable for wall-mounted lasers.

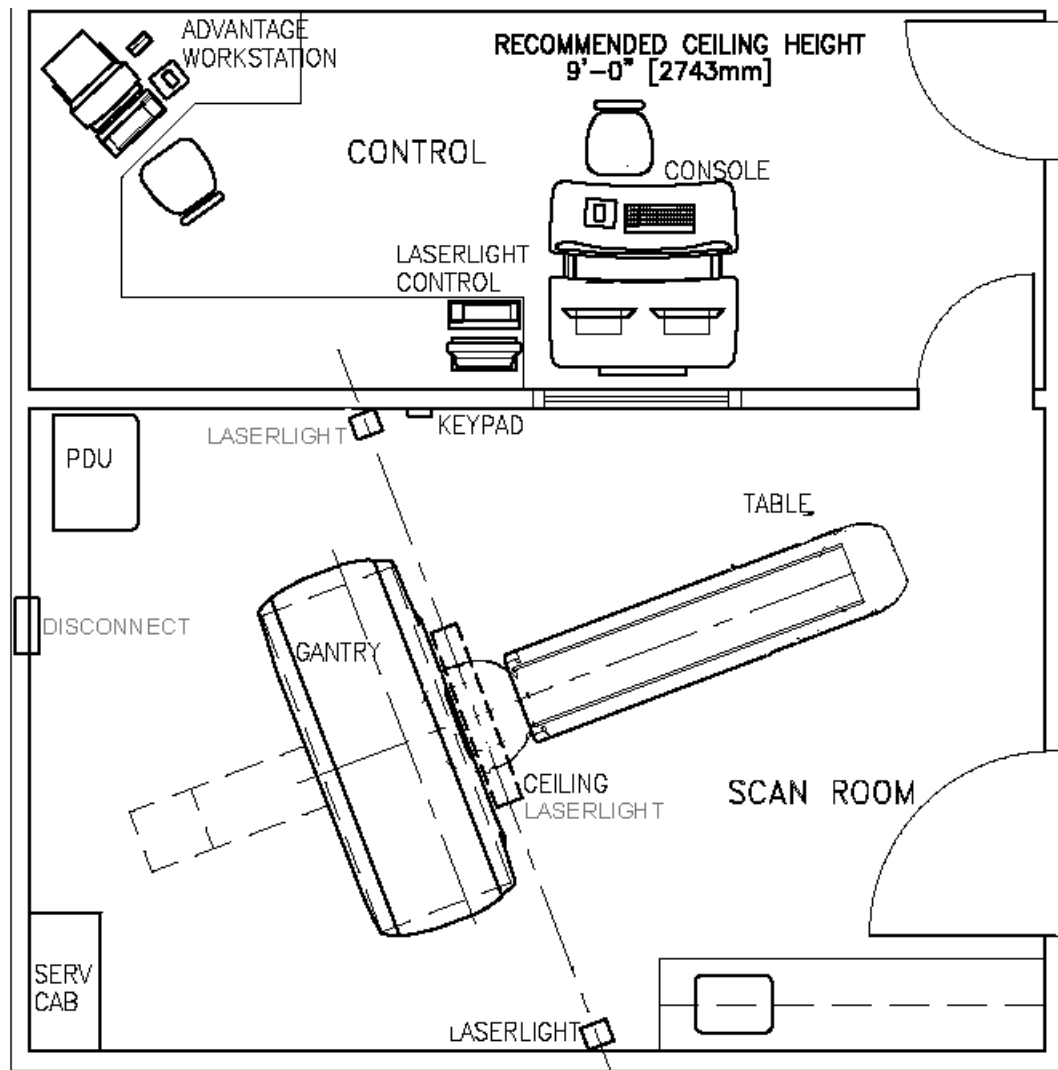


Figure 4-6 Typical Room Layout with CT Simulation Laser Alignment System

Section 6.0

Component Dimensions

DESCRIPTION	WIDTH		DEPTH		HEIGHT	
	MM	INCH	MM	INCH	MM	INCH
Gantry	2439	96.0	1006	39.6	1992	78.4
Table	610	24	2387	94		
Power Distribution Unit (NGPDU)	711	28	559	22	1067	42
Operator's console/computer	1219	48	991	39	851	33.5
Remote Color Monitor (LCD)	413	16.25	203	8	406	16
Color printer	584	23	457	18	178	7
Workstation	419	16.5	409	16.1	79	3.1
Monitor	495	19.5	483	19	471	18.6

Table 4-8 Dimensions of Components

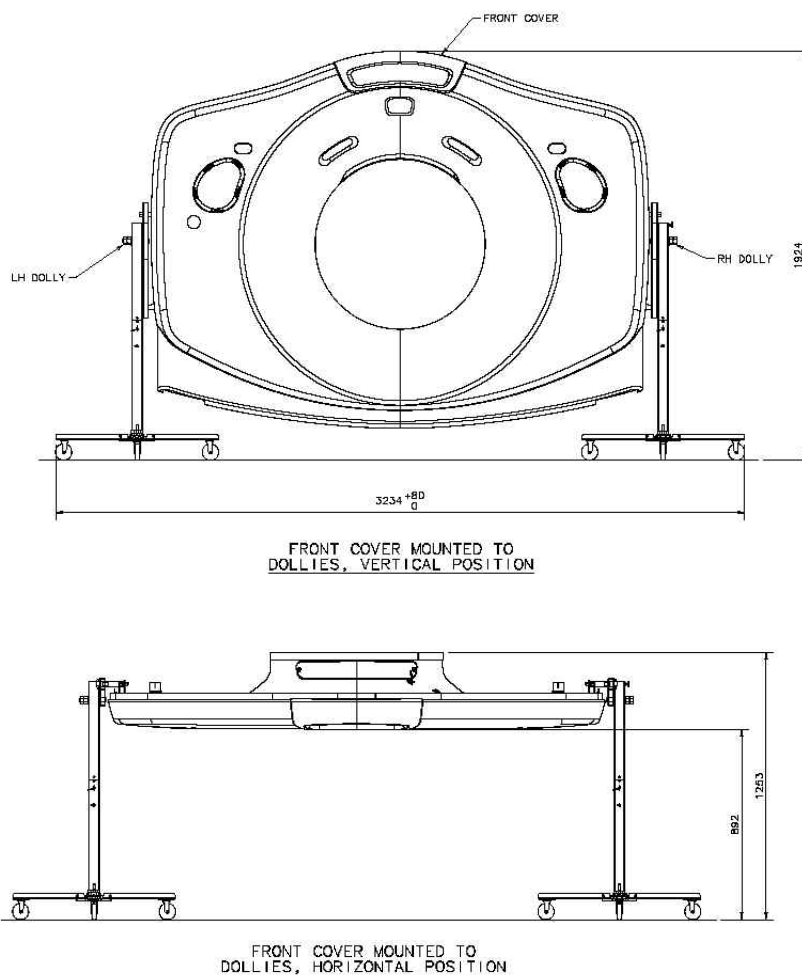


Figure 4-7 Gantry Front Cover with Service Dolly Dimensions

6.1 Table and Gantry

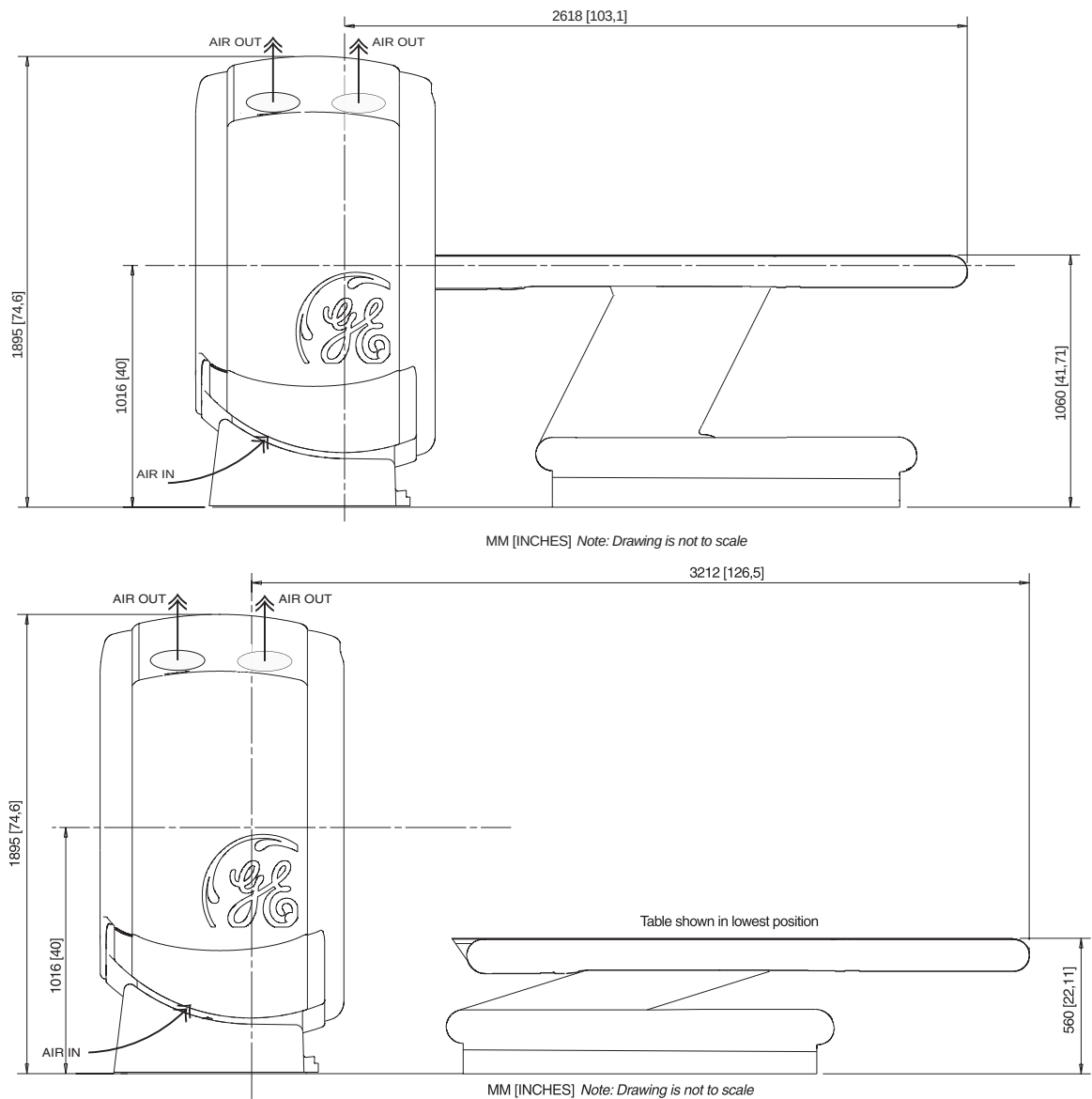


Figure 4-8 Table and Gantry (Side View)

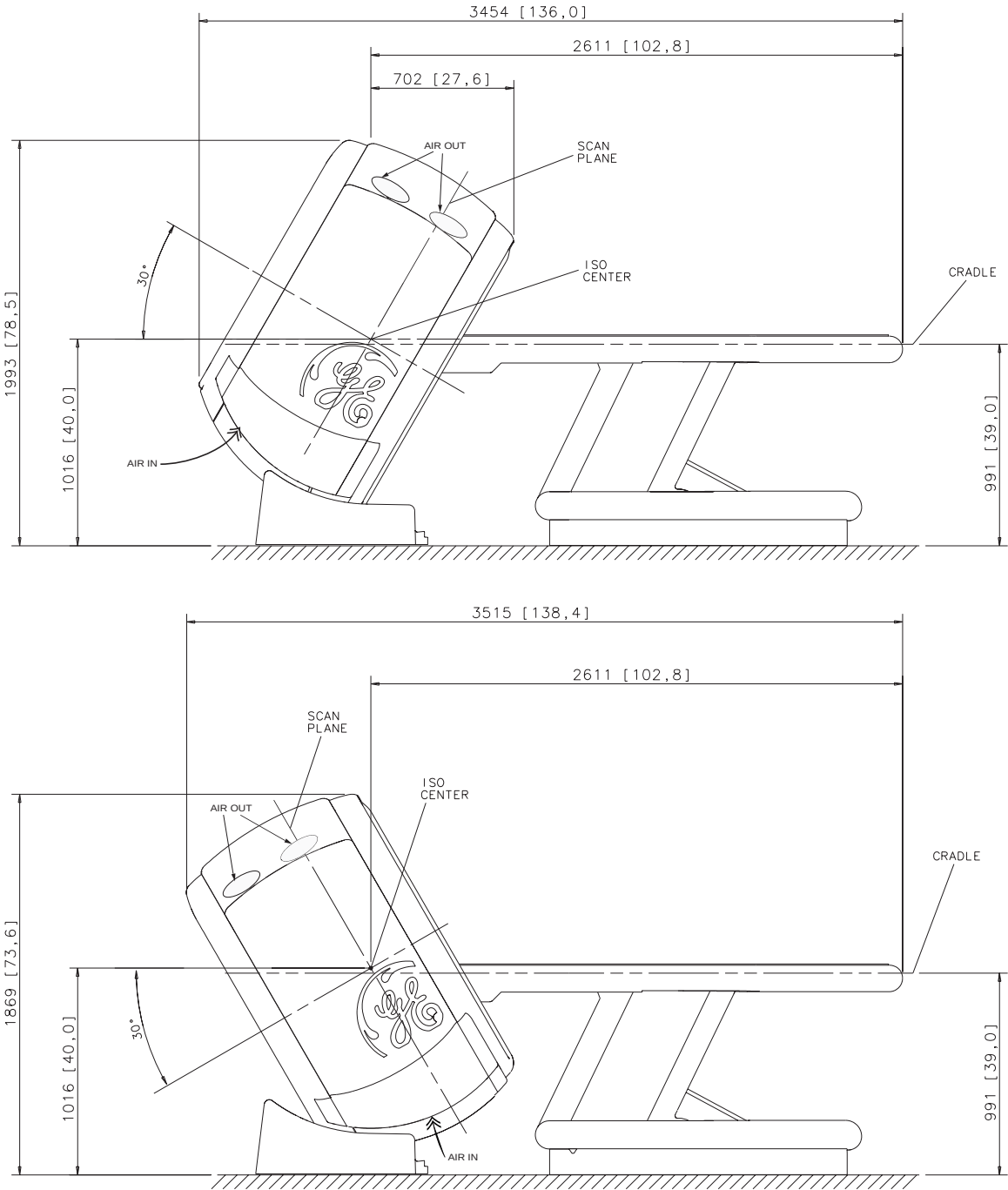


Figure 4-9 Gantry shown tilted +30° (top) and -30° (bottom)

6.2 Power Distribution Unit

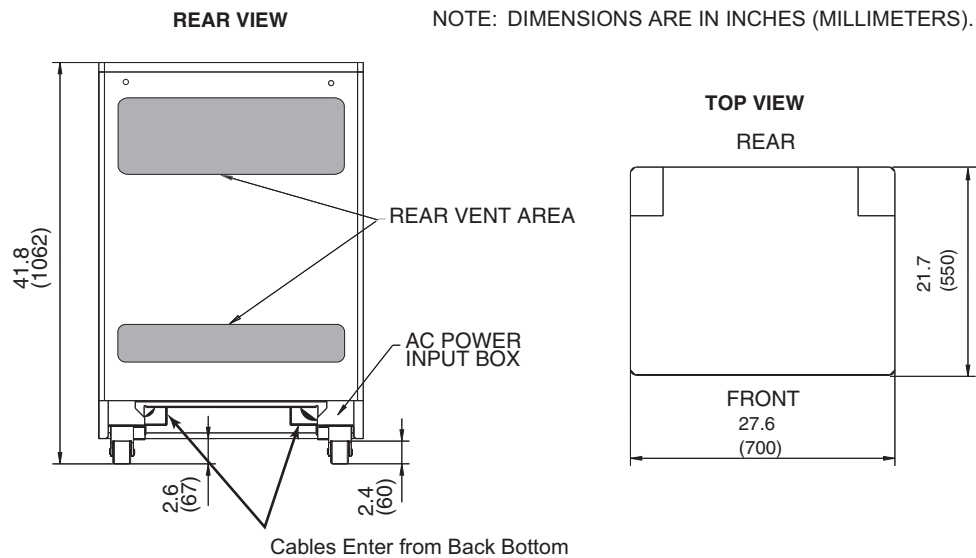


Figure 4-10 Power Distribution Unit (NGPDU-3)

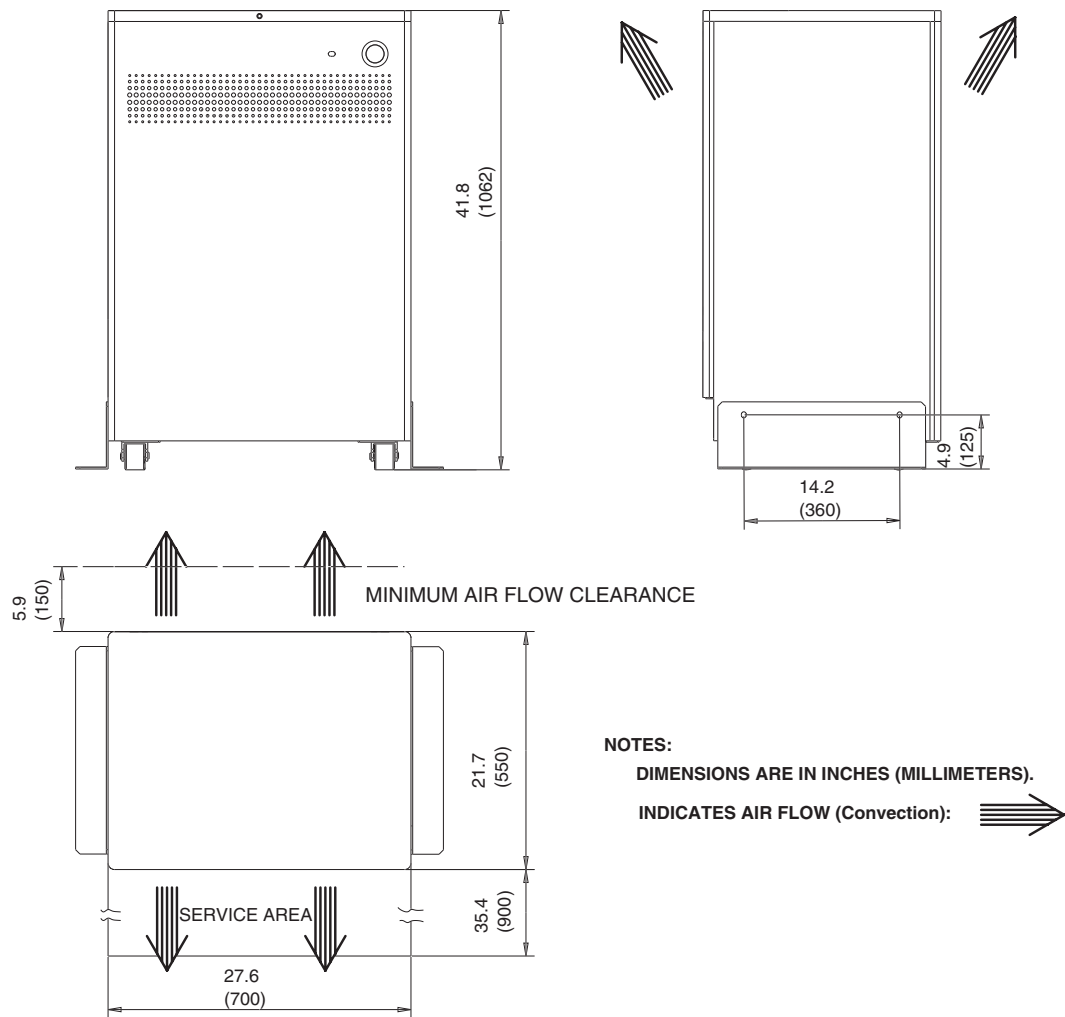


Figure 4-11 Power Distribution Unit (NGPDU-4)

6.3 Operator's Console

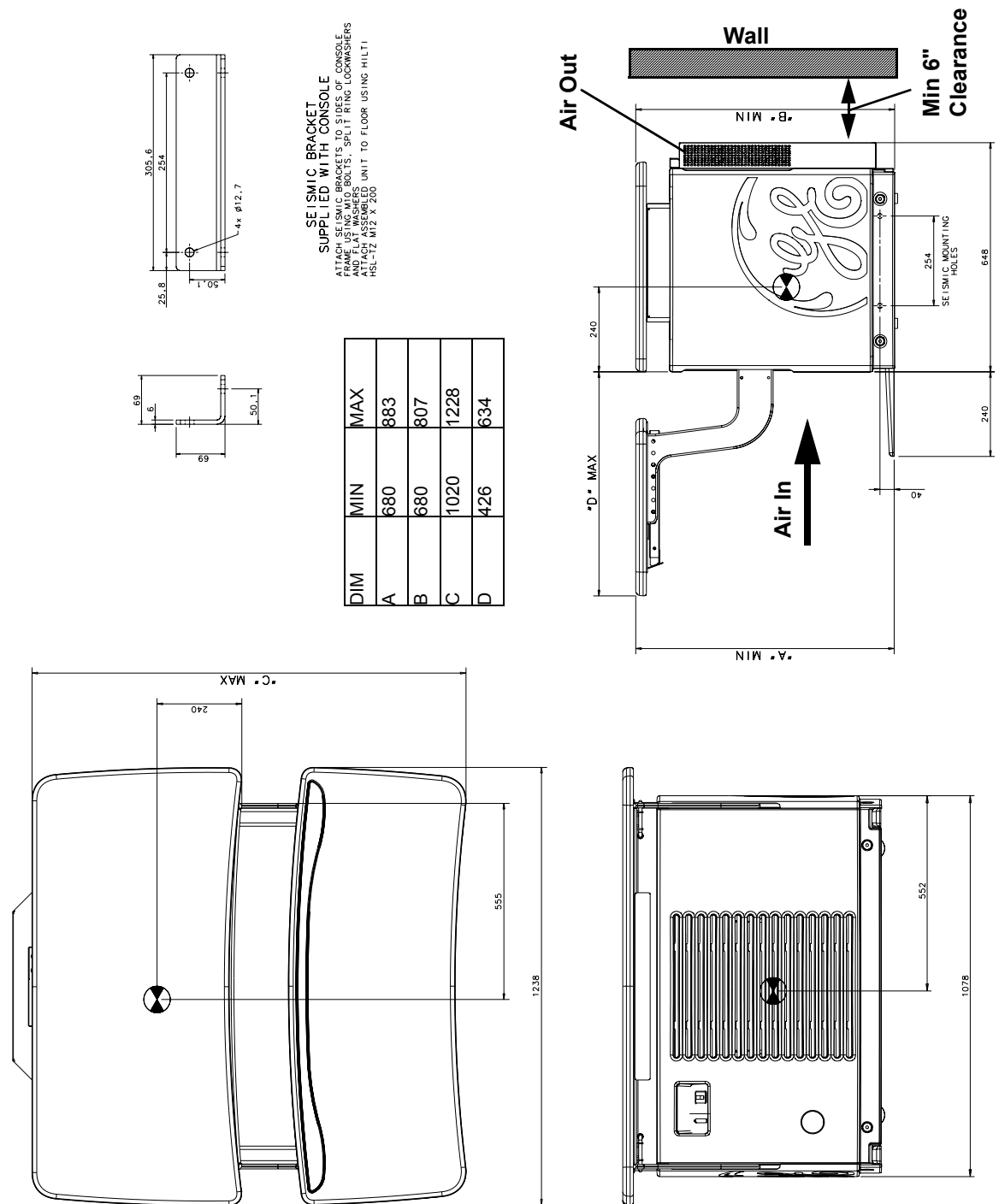


Figure 4-12 Operator's Console

Section 7.0

Minimum Dimensions and Clearances

SYSTEM OPERATION	MM	INCHES
Ceiling Pedestal mount (optional) (lowest point to floor injector or monitor)	2438.5 mm	96 in.
Finished ceiling to floor (recommended)	2743 mm	108 in.
Finished ceiling to floor (minimum)	2438.5 mm	96 in.
Table max extension head end w/ extender from Center Line	2030 mm	80 in.
Table extension head end w/extender to obstruction	152 mm	6 in.
Table in lowest pos. w/cradle @home pos. to Center Line	3209 mm	126.5 in.
Back of Console to wall	152 mm	6 in.
Back of PDU to wall	152 mm	6 in.

Table 4-9 Minimum Dimensions and Clearances

7.1 Injector Control

A suitable work area, which is within reach of the Operator's Console, should be provided for placement of the injector control. Refer to [Figure 4-13](#).

Wall mounted, ceiling mounted and pedestal units need cables to be routed from the Gantry area to the Operator Console area. The supplied cable is 15.24 m (50 ft.) in length. Injectors require AC power that is not supplied by the CT system.

Note: Injector cables should not be routed with the system cables.

Mounts are available in different configurations and lengths. Refer to Injector documentation for detailed installation instructions.

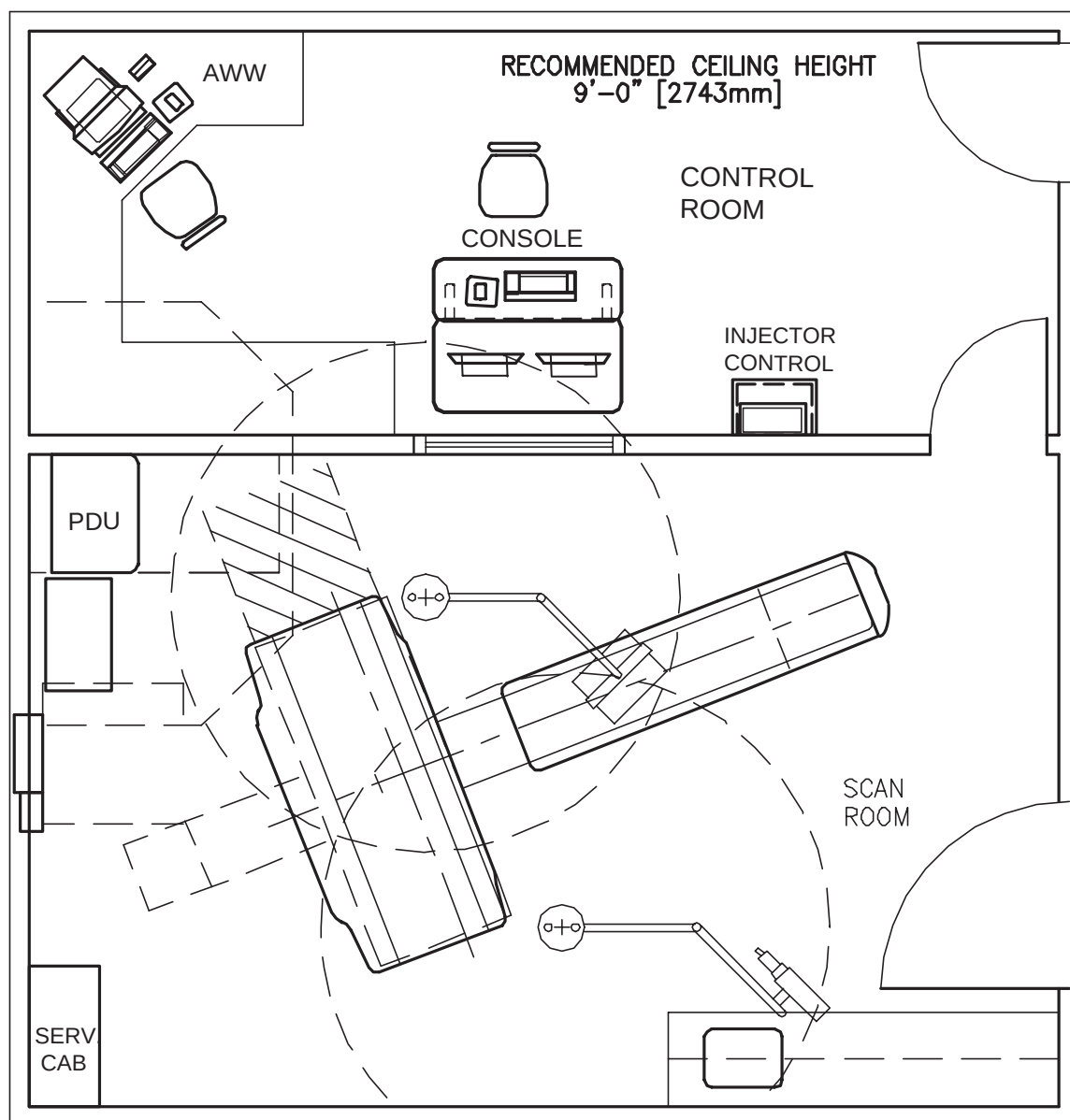


Figure 4-13 Typical Room Layout, with Injector Option

Section 8.0

Structural Requirements

8.1 Table and Gantry Mounting Requirements



WARNING **POTENTIAL FOR PATIENT INJURY.**
IMPROPERLY SECURED TABLE MAY TIP, DISLODGING PATIENT.
PROPER ANCHORING IS KEY TO MAINTAINING PATIENT SAFETY DURING SYSTEM OPERATION.



NOTICE **It is the purchaser's responsibility to provide an approved support structure and mounting method for all floor types other than those listed. General Electric is not responsible for any failure of the support structure or method of anchoring, including seismic requirements.**
GE is not responsible for methods other than those listed.

Table and gantry mounting dimensions are shown in [Figure 6-2](#), [Figure 6-3](#), [Figure 6-4](#). Refer to [Chapter 6](#) for additional details of floor loadings, component weights, and Gantry and Table installation and anchoring.

Anchor gantry and table to floor by a means that will maintain their relative alignment and meet applicable building and other local codes, including seismic structural mounting requirements.

Floor structure must be capable of withstanding the occupied weight of table and gantry, and the individual contact area loading of these components.

Support areas of the patient table and gantry must rest on solid concrete, not resilient tile or carpeting which will slowly yield over a period of time and disturb alignment of table to gantry.

Factors that could cause misalignment between gantry and table due to floor sag should be considered. The cradle can potentially carry a 204kg (450 lb) patient. Center of gravity changes as cradle cantilevers.

Take into consideration all other moving weights such as gurneys or personal equipment. Refer to [Chapter 6](#) for gantry and table mounting details.

No part of floor surface within table and gantry, nor the two interface areas between table and gantry, should be higher than the support area for table and gantry.

8.2 Floor Anchors

Provided floor anchors are designed for use ONLY on concrete floors that meet the 101.6 mm (4 in.) concrete floor requirements. All other anchoring methods (on floor types other than the 101.6 mm (4 in.) concrete minimum) must be determined to meet the stated GE minimum load requirements, at the customer's expense, by their structural contractor. The customer's contractor is responsible for the installation of all anchors other than those shipped with the system

8.3 Floor Strength

Concrete floors must have a minimum strength of $f'_c = 2000$ psi (1.4×10^7 Pa) at 28 days for (curing time) mounting floor anchors. It is the responsibility of each customer to have appropriate tests performed to determine and measure concrete strength.

Note: If installing the GE LS scanner on a floor type thinner than a 101.6 mm (4in.) concrete floor, the purchaser, at their expense, shall provide acceptable anchoring and mounting methods that meet all structural specifications provided in sections 8.1 through 8.5 of this Pre-Installation Manual.

8.4 Floor Levelness

The CT Room floor levelness requirement is important for accurate patient positioning. Floor levelness in the Scan Room must not be greater than 3 mm (0.125 in.) between depression and high spots over any 3048 mm (120 in.) distance within the area of the gantry/table template (see the envelope shown in [Figure 2-1, on page 43](#)).

Note: The floor must meet levelness specification to properly align the table/gantry. Minimum gantry height at this specification is 15 mm (1/2 in.) to prevent cable crushing.

Table level may not be achievable if the overall floor levelness is greater than the specification. The overall floor level must be 0 in. to use under gantry cable entrance. The minimum gantry height is 20 mm (3/4 in.) with this option to prevent cable crushing.

The use of floor shims is not suitable to achieve floor levelness. It is recommended that the concrete be leveled to meet this requirement.

8.5 Floor Vibration

The CT equipment may be sensitive to vibration in the frequency range of 0.5 to 20 Hz depending on the amplitude of the vibration. It is the customers responsibility to contract a vibration consultant or qualified engineer to implement design modifications to meet the specific limits, However, it is ultimately the customer/architect/engineer responsibility to design the site solution.

8.5.1 Steady State Vibration

The maximum steady state vibration transmitted through the floor should not exceed 10^{-3} m/s² rms maximum single frequency above ambient baseline from 0.5 to 80 Hz (measured in any 1 hour during a normal operating period).

8.5.2 Transient Vibration

The behavioral characteristics must be such that any measurable transient disturbance must also be minimized to less than 0.01 m/s² peak-to-peak.

8.5.3 Equipment Location

To minimize the interference, the CT equipment should be placed on a solid floor, located as far as possible from the following vibration sources:

- Parking lots
- Roadways
- Subways
- Trains
- Hallways
- Elevators
- Heliports
- Hospital power plants containing pumps, motors, air handling equipment and air conditioning units

8.6 Walls: Scan Window

The recommended patient viewing window dimensions are 1219 mm wide x 1067 mm high (48 in. x 42 in.). The location of the window is dependent on the position of Operator Workspace position. Consult [Section 10.0](#) of this chapter and a **qualified radiological health physicist** for radiation protection requirements of glass (lead content and thickness).

Note: The operator at the Operator Workspace must be able to view the patient during a scan.

Section 9.0 Network Connections

Broad-band is the standard network connection and is required for all LightSpeed 5.X systems. (A dial-up modem is optional.) A 1000 baseT high-speed network is desired, with 100 baseT network service acceptable. Broad-band connections should use one of the following Category 5 patch cables:

CAT Num	GE Part Num	Length
K9000WB	2215028-10	20 m
K9000KP	2215028-5	10 m
K9000JR	2215028-4	5 m
K9000WA	2215028-9	3 m

The CT system is connected to the network through the Console.

- A patch cable (not to exceed 3.05 m (10 ft.)) should be provided by the customer, and it is used to connect the console to a wall box. (See Notes on [Figure 9-2](#))
- Some customer-site units may require cable duct-work or conduit to route connecting network cables to the workstation, camera and console.
- The run from the hospital switch to the CT wall outlet must not exceed 88 m (290 ft.). Bandwidth performance is degraded when the length reaches 91 m (300 ft.) or greater.
- For the optional modem: **Two phone lines should be provided by the facility.** One line is for

use with a modem and must be an analog line. The second line is a voice only line.



Figure 4-14 Console Rear Bulkhead

9.1 US Broad Band Process Overview

The United States network connectivity requirement for this product is broad-band. The US process relies on the Project Manager Installation to select a Customer Champion and identify an IT contact for the site. Together, those individuals then complete a site assessment to gauge what tasks are needed to fulfill the connection.

Anyone can contact the GE Connectivity team at 800.321.7937, Option #3, with questions.

9.2 Customer Broad-Band Responsibilities

Provide GE Healthcare Project Manager of Installation (PMI) with an accurate site address, telephone number, contact name, and e-mail address for the:

- Customer Champion
 - Coordinate VPN activities between Radiology/Cardiology and the Information Technology (IT) departments
 - Act as a focal point in assuring site broad-band infrastructure meets GE Healthcare requirements for connection as determined by a mutual assessment with the GE Healthcare Connectivity team.
- IT Contact
 - Complete an equipment assessment with GE Healthcare Connectivity team to determine site readiness for broad-band
 - Contact your Project Manager of Installation (PMI), for the name of the zone Broad-Band Specialist
 - Work with the Customer Champion to complete any identified infrastructure changes
 - Provide IP addresses for new CT equipment
 - Provide a VPN compatible appliance that will support the IPSec tunneling protocol and 3DES data encryption
 - To utilize an Internet Service Provider that supports static routing



GE Medical Systems

Remote Service Broadband - Customer Site Assessment

Site Name: _____ FE Name: _____
City, State: _____ FE Phone: _____
Date: _____ FE Email: _____

- | | Yes | No |
|--|--------------------------|--------------------------|
| 1. Does your site currently have a persistent (24x7) Internet connection? | ☐ | ☐ |
| 2. Is the GEMS Diagnostic Imaging equipment on the Local Area Network and will it be accessible to the Internet? | ☐ | ☐ |
| 3. Does your site have a VPN device today? | ☐ | ☐ |
| 4. Is the VPN device one of the models below? <i>If Yes, please select the model from the options below.</i> | ☐ | ☐ |

- | | |
|--|---|
| ☐ a) Cisco Pix Firewalls | ☐ g) Symantec (Raptor) firewalls |
| ☐ b) Cisco Routers | ☐ h) Firebox |
| ☐ c) Cisco 3000 Series (Altiga acquisition) | ☐ i) Linux S/WAN |
| ☐ d) Checkpoint Firewalls Software Version 4.1 and higher | ☐ j) Sidewinder |
| ☐ e) Nortel Contivity Software Version 3.2 or higher | ☐ k) Netscreen |
| ☐ f) Redcreek | ☐ l) None |
| | ☐ m) Other _____ |

**If No, the GEMS Connectivity Support Team can help determine device compatibility.*

- | | | |
|--|--------------------------|--------------------------|
| 5. Does your VPN device support "triple DES" Encryption? | ☐ | ☐ |
| 6. Has approval been given to install this VPN connection? | ☐ | ☐ |
| Site Approver's Name _____ | | |
| 7. Provide your VPN Installer information, this is the person who will be contacted to schedule the VPN install. | | |

Customer Installer Name: _____
Installer Telephone Number: _____
Installer e-mail address: _____

Notes: _____

Field Engineer needs to provide compatible system information:

System ID	IP Address	Gateway Address
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____
_____	_____	_____

[Additional System and IP Address Spaces Available on Page 2](#)

If you have questions or need assessment support contact your Zone Champ or:

Joe Gracz - HQ Support 1-262-524-5261
Joseph.Gracz@med.ge.com

Once you have completed both pages of this form, please send it to:

a) Judy Heyer judy.hey@med.ge.com
b) Judy Heyer Fax# 414-918-4707

[Use the send button on page 2](#)



System ID List & Compatibility Matrix


MR	Software Version	CT	Software Version
Include all Signa(LX, Ovation, OpenSpeed and CV/I) MR's System Ids and IP Addresses in assessment.		HiSpeed CT/e	Minimum 6.0
		HiSpeed CT/e dual	Minimum 2.06
		HiSpeedNX/I	Minimum 5.5
		LightSpeed QX/I	Minimum 1.3
		LightSpeed Plus	Minimum 2.1
Prior to completing an InSite checkout systems must be upgraded to Version 9.0 or 9.1 unless shown in the compatible list below.		LightSpeed Ultra	Minimum 3
Signa TwinSpeed	Release 9.0	HiSpeed X/I	
Signa	Release 9.1	CT/I	Minimum 6.2
Signa	Release 10.x		
Signa CV/I	CNV4	Lxi, Dxi, Fxi, Zxi	Minimum 6.01
Signa 3T			
Signa Profile	Minimum 7.66		
Signa Contour	Minimum 7.66		
Magnet Monitor	Minimum 2.3 software		
Network Products	Software Version	Nuclear/Pet	Software Version
AW 4.0	Version 4.0 or Above	Advance	Minimum 5.1
AW 4.0P		eNTEGRA	Minimum 2.03
X-Ray	Software Version		
<u>Cardiac</u>			
INNOVA	All		
<u>Mammography</u>			
Seno2000D	All		
<u>Digital Radiology</u>			
Revolution XQ/I	Version 10.12.5 or above		
Revolution XR/d	Version 18.0		
Multi Vendor			
MR	PSI	CT	PSI
Philips ACS 1.5T	MZP 400	Picker MX-TWIN	CZM 400
Philips NT 2000	MZP 401	Picker MX 8000	CZM 401
Philips NT 3000	MZP 402	Picker PQ	CZM 500
Philips NT ACS	MZP 403	Picker PQ 2000	CZM 501
Picker Outlook 0.23T	MZM 301	Picker PQ 5000	CZM 502
Picker Edge 1.5T	MZM 800	Picker PQ 6000	CZM 503
Picker Vista 1.0T	MZM 801	Picker PQS	CZM 504
Picker Eclipse 1.5T	MZM 900	Siemens AR STAR	CZS 100
Picker Polaris 1.0T	MZM 901	Siemens AR.C	CZS 101
Siemens Harmony 1.0T	MZS 100	Siemens AR.HP	CZS 102
Siemens Impact 1.0T	MZS 200	Siemens AR.SP	CZS 103
Siemens Impact Expert 1.0T	MZS 201	Siemens AR.T	CZS 104
Siemens SP 1.0T	MZS 400	Siemens PLUS-4 EXP	CZS 400
Siemens SP 1.5T	MZS 401	Siemens PLUS-4 POW	CZS 401
Siemens SP4000 1.5T	MZS 403	Siemens PLUS-4 VOL	CZS 402
Siemens Symphony 1.5T	MZS 500		
Siemens Vision 1.5T	MZS 600		

Field Engineer needs to provide compatible system information:

[illegible]

Once you have completed this form, please send it to:

judy.heyer@med.ge.com
414-918-4707

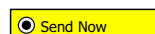

Section 10.0 Radiation Protection



NOTICE Scanner-room shielding requirements should be reviewed by a qualified radiological health physicist taking into consideration:

- Scatter radiation levels within the scanning room (see [Figure 4-15](#))
- Equipment placement
- Weekly projected work-loads (# patient/day technique (kvp*ma))
- Materials used for construction of walls, floors, ceiling, doors, and windows
- Access to surrounding scan room areas
- Equipment in surrounding scan room areas (e.g., film developer, film storage)

[Figure 4-15](#) depicts measurable radiation levels within the scanning room while scanning a 32 cm CTDI phantom (body) and a 20 cm water phantom (head) with the technique shown.

Note: Actual measurements can vary. All measurements have an accuracy of $\pm 20\%$ because of measurement equipment, technique, and system-to-system variation.

Use the correction factors shown in [Table 4-10](#) to adjust exposure levels to the usual scan technique at your site.

CHANGED PARAMETER	MULTIPLICATION FACTOR
mAs	new mAs/100
80 kV	0.21
120 kV	0.71
140 kV	1.0
4 x 3.75mm images	0.82
16 x 0.625 LD	0.59
8 x 1.25 LD	
Fluro 5mm	
4 x 1.25 LD	0.40
5mm (1i)	
Fluro 2.5 mm	
1 x 1.25mm images	0.20
4 x 0.625 LD	0.10
1 x 1.25	

Table 4-10 Shielding Requirements Scaling



NOTICE The units of measure used for radiation levels have been changed in this publication, from mR (millirads) to μGy (micrograys). The conversion factor is:

$$1 \text{ mR} = 8.69 \mu\text{Gy}$$

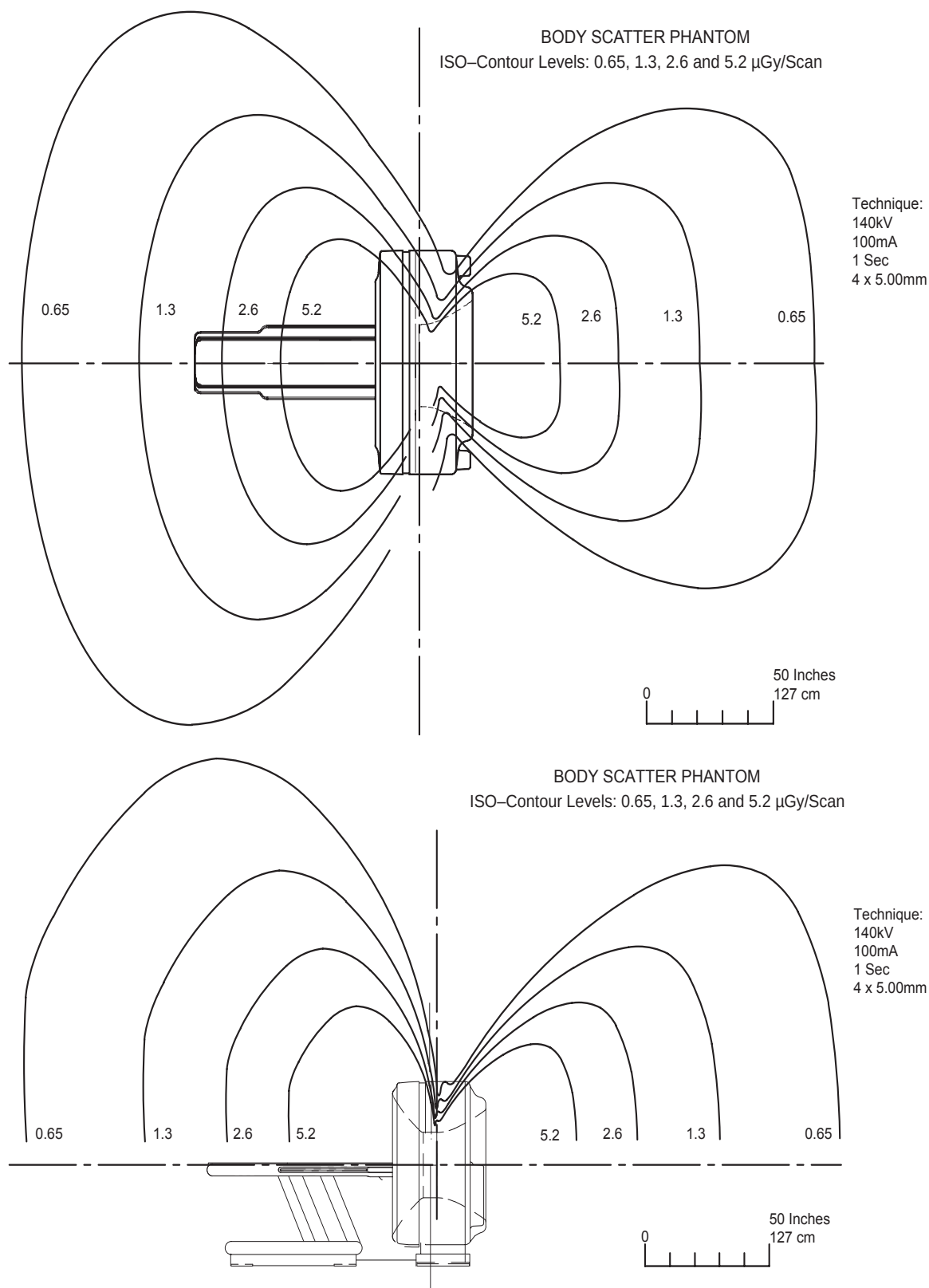


Figure 4-15 Typical Scatter Survey (Body Phantom)

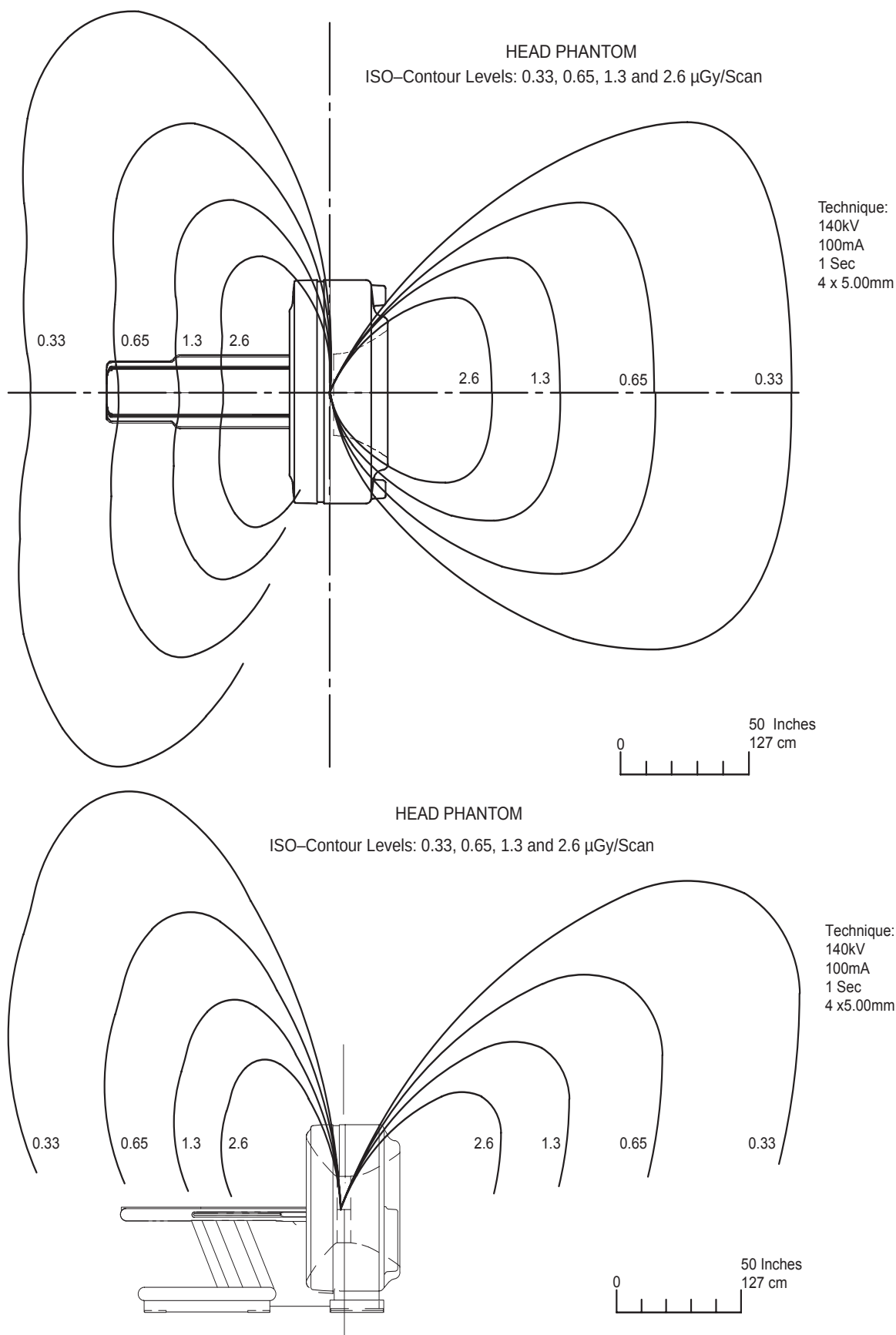


Figure 4-16 Typical Scatter Survey (Head Phantom)

Chapter 5

Environmental Conditions

Ratings and duty cycles of CT subsystems apply if site environment meets the standards of this section. Maintain environmental conditions listed below at all times – including, for example, overnight, weekends and holidays. Shut down the CT system if air conditioning is not working. When system is shut down for major repair, air conditioning may be shut down also.

Section 1.0

Temperature and Humidity Specifications

Environmental specifications apply to the table, gantry, power distribution unit, and operating console.



NOTICE System Operation and Image Quality may be affected if environmental specifications are exceeded.

1.1 Temperature (Scan Room & Control Room)

Maximum allowable ambient room temperature:	26° C (79° F)
Recommended ambient room temperature:	22° C (72° F)
Minimum allowable ambient room temperature:	18° C (64° F)
Maximum allowable ambient room temperature rate of change: 3° C per hour (5° F per hour)	

Note: Any cooling equipment cycle control range must be taken into account, such that the maximum and minimum ambient room temperatures shown above are not exceeded, during room thermal cycling. For example, if the HVAC is capable of ± 2° C control, then the limits would 20° C - 24° C, to maintain absolute limits.

1.2 Humidity (Scan Room & Control Room)

Maximum allowable non-condensing relative humidity:	60%
Minimum allowable non-condensing relative humidity:	30%
Maximum allowable relative humidity rate of change:	5% per hour

1.3 Other Guidelines

- To help determine the hospital room environmental conditions, a temperature and humidity recorder may be temporarily installed close to where the system will be installed. Note readings before and after installation, to verify the true temperature and humidity in your environment.
- Consider heating, ventilating, air conditioning (HVAC) needs and redundancy. An air conditioner with two compressor units, rather than one, may prevent system downtime. A back-up (redundant) air conditioner permits CT system operation during an extended repair of the primary air conditioner.

Section 2.0 Cooling Requirement

Use [Table 5-1](#) to assist in cooling requirements planning. Over half the cooling used by your scanner is required for gantry operation. Locate a wall air-conditioning vent at floor level beside and behind gantry to meet both gantry cooling needs and provide patient comfort. Do not locate any cooling vents directly above the gantry. Air returns above the gantry are recommended.

SYSTEM COMPONENT	BTU/HR	WATT
1. Gantry recommended (See NOTE 1)	43,000	12,600
2. Table	700	200
3. Power distribution unit	3400	1000
Recommended Scan Room Subtotal (see notes):	47,100	13,800
4. Operator's console/computer with one IG	7361	2165
Additional IG, each	1360	400
LCD Monitor, each (x2)	170 (340)	50 (100)
SCSI Tower	425	125
Injector	425	125
DASM	425	125
Recommended Control Room Subtotal (includes OC w/1 IG, 2 monitors & SCSI Tower):	8126	2390
System Total (Recommended) (See NOTE 1)	55,226	16,190
Option: Remote Color Monitor (LCD)	170	50
Option: UPS	2900	850
ROOM TOTAL (SEE NOTE 2)		
NOTE 1: With 75 scan rotations per patient: Recommended BTU/hr. provides for up to six patients per hour. It is also needed during calibration of the system.		
NOTE 2: Cooling requirements do not include cooling for room lighting, personnel or non-CT equipment.		

Table 5-1 Cooling Requirements (Worksheet)

Section 3.0

Altitude

System operating altitude is from mean sea level to 10,000 ft. (3050 meters).

Section 4.0

Electro-Magnetic Interference (EMI)

Note: If power sub-stations exist under or above the scan room, or near the control room, consider EMI testing to determine if your proposed room meets the published acceptable EMI room limits. This also includes high voltage lines under the scan or control room floor.

4.1 Gantry

Locate gantry in ambient static magnetic fields of less than 10^{-4} tesla (1,000 milligauss) to guarantee specified imaging performance. Ambient AC magnetic fields must be below 10^{-6} tesla (10 milligauss) peak.

4.2 Color Monitor

Locate color monitors in ambient static magnetic fields of less than 5×10^{-5} tesla (100 milligauss) to guarantee color purity and display geometry. See [Figure 5-1](#).

4.3 Console / Computer Equipment

Locate computer equipment in ambient static magnetic fields of less than 10^{-3} tesla (10,000 milligauss) to guarantee data integrity. See [Figure 5-1](#).

4.4 Magnetic Media

Locate magnetic media in ambient static magnetic fields of less than 10^{-3} tesla (10,000 milligauss).

4.5 PDU

The PDU produces an electromagnetic field that radiates outward from its cabinet in all directions. Do not place sensitive electronics (e.g., console or computer equipment - the UPS is not classified as sensitive electronics) within one-and-a-half meters (1.5m) of the Power Distribution Unit., in any direction (including above or below). See [Figure 5-1](#).

4.6 EMI Reduction

If fields of excessive EMI are known or suspected to be present, consult GE Healthcare Sales & Service for recommendations. Consider the following if you attempt to reduce EMI:

- External field strength decreases rapidly with distance from source of magnetic field.
- External leakage magnetic field of a three-phase transformer is much less than that of a bank

of three single phase transformers of equivalent power rating.

- Large electric motors are a source of substantial EMI.
- High-powered radio signals are a source of EMI.

Maintain good screening of cables and cabinets.

4.7 UPS

The Uninterruptable Power Supply (UPS) provides a consistent power supply to various electrical components of the system. Also, it continues to provide electrical power to components during a site-wide power outage so components can be safely shut down. The UPS should be kept at least one meter (1m) away from sensitive electronics (the PDU does not include sensitive electronics).

For UPS interconnect information, please refer to [Chapter 9, Section 6.0](#) on [page 116](#).

4.8 Equipment EMI “Envelopes”

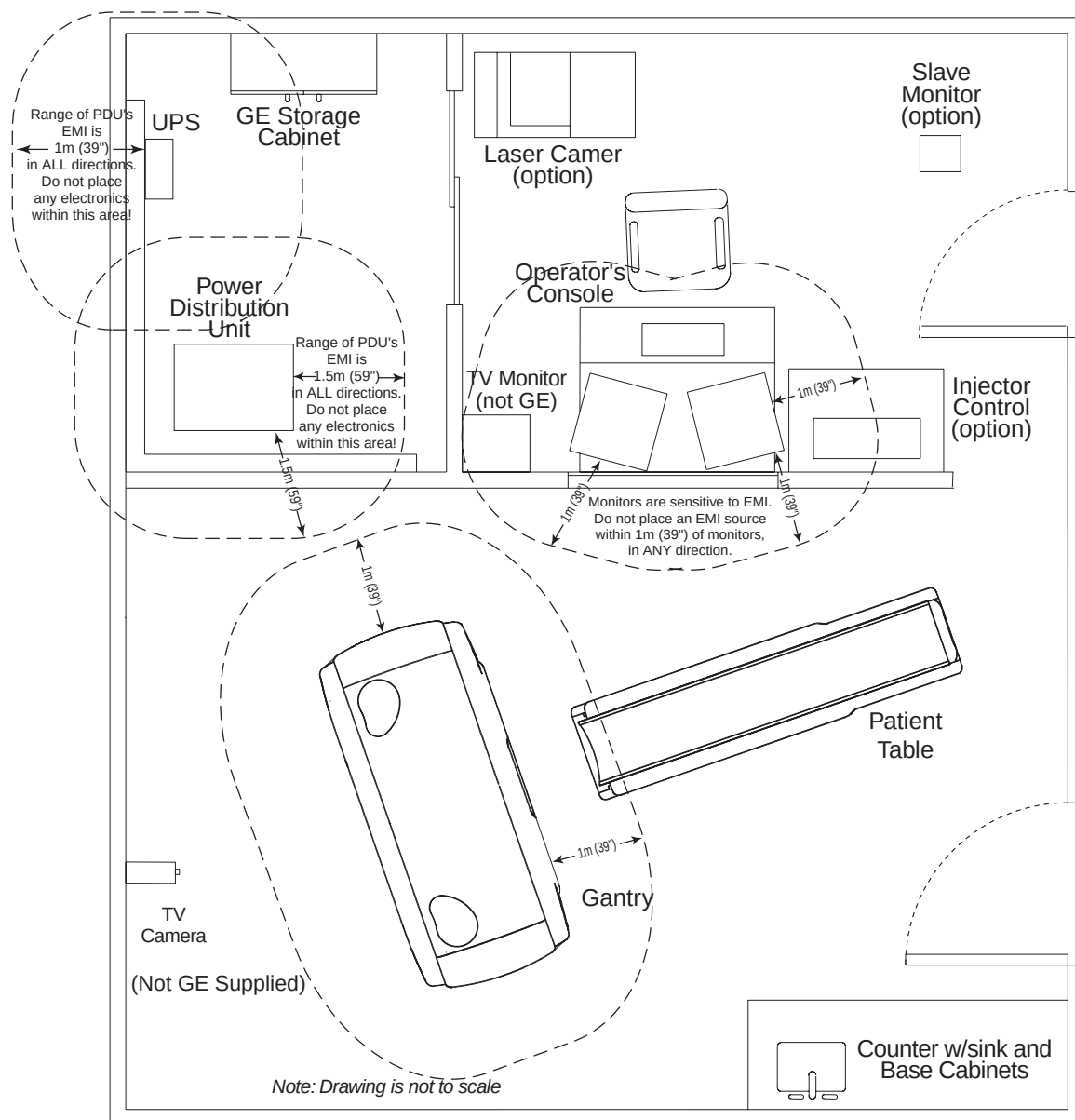


Figure 5-1 Sample Room Layout, showing approximate EMI requirements

Chapter 6

Floor Loading and Weights

Section 1.0 Floor Loading

The LightSpeed RT system has a total floor load of approximately 6660 lbs (3021 kg). About 5485 lbs (2500 kg), including patient (450 lbs (204 kg)), is concentrated in the table-gantry assembly.

Table 6-1 lists CT components with weight, size, floor loading and normal mounting requirements.

ITEM	NET WEIGHT LB(KG)	OVERALL W X D INCH (MM)	WEIGHT/AREA LB/SQ. FT. (KG/M ²)	LOAD PATTERN IN. (MM)	NORMAL METHOD OF MOUNTING IN. (MM) (GE SUPPLIED) ¹
Gantry (w/covers)	~3891 (~1765)	96.0 X 39.6 (2439 X 1006)	238 (1164)	Rectangular base plate 24 X 81 (610 X 2057) with four round pads, each 2.5 (63.5) in contact with floor.	Hilti Kwik-Bolt II 1/2in (12.7mm) diameter by 8in (203mm) long per P/N 2106573 at four leveling pads into concrete floor.
Dollies (each)	250 (114)			Individual pad loadings are: 910 lb, 960 lb, 1040 lb, and 1090 lb (see Figure 6-3).	
Top Cover (each)	24.6 (11.2)				
Side Cover (each)	25 (11.3)				
Front Cover	95 (43)				
Rear Cover	100 (45)				
Patient Table	1185 (538) Includes 450 (204) Patient	24 X 94 (610 X 2387)	175 (852)	Rectangular base 17 X 57.5 (432 X 1460) with five round pads, each 2.5 (63.5) in contact with floor. Worst-case cantilever pull on any bolt is 1235 lb (561 kg).	Hilti Kwik-Bolt II 1/2in (12.7mm) diameter by 8in (203mm) long per P/N 2106573 at three leveling pads into concrete floor.
Dolly	215 (98)				
Power Distribution Unit	~700 (~336)	28 X 22 (711 X 559)	161 (788)	Four Casters support area of 28 X 22 (711 X 559).	Casters are for positioning and service. Set on floor. May be anchored to floor using angle brackets ² in seismic zones.
Console w/HP & w/o monitors	450 (204)	48 X 39 (1219x991)	140 (681)	Four Casters or Leveling Feet support area of 46 X 19 (1168 X 483).	Casters are for positioning. Set on floor. Console may be anchored to floor using angle brackets ²
Monitor - CRT (ea)	80 (36)				
Monitor - LCD (ea)	22 (10)	16.3 X 17.8 (415 X 453)			

Notes:

- 1.) Use the GE Supplied mounting hardware ONLY IF APPROVED by qualified personnel.
[See statements in 8.1 - Table and Gantry Mounting Requirements, on page 72.]
- 2.) Angle brackets are included on the shipping skid, and are also available in the OPTIONAL Seismic Kit (R4390JC).

Table 6-1 LightSpeed RT System Floor Loading

Section 2.0

Mounting Data, Including Seismic

 **WARNING POTENTIAL FOR PATIENT INJURY.**
IMPROPERLY SECURED TABLE MAY TIP, DISLODGING PATIENT.
PROPER ANCHORING IS KEY TO MAINTAINING PATIENT SAFETY DURING SYSTEM OPERATION.

The following pages show center-of-gravity information for system components:

- Gantry: [Figure 6-2](#)
- Table: [Figure 6-4](#)
- Power Distribution Unit: [Figure 6-5](#)
- Operator's Console/Computer: [Figure 6-7](#)

Floor mounting hole locations for components that don't have templates are also in this section.

An optional Seismic Kit (R4390JC) is available, for facilities in areas that require one.

Customer is responsible for seismic mounting. Refer to all applicable codes in your area.

GE provided floor anchors ([Figure 6-1](#)) are designed to be used ONLY on concrete floors that meet the concrete floor requirement. Supplied floor anchors must be installed by a trained contractor, and shall be set to a minimum depth of 3 inches at each anchor point. Any anchors having more than 1 inch of thread showing above the nut, when torqued to 55 ft.-lbs, shall have a second anchor installed in the closest adjacent hole. the second anchor shall be installed to the standard depth and torque specifications.

MOUNTING REQUIREMENTS	GANTRY	TABLE
Minimum Floor Thickness:	4 inches	4 inches
Recommended Drilling Depth:	3-¾ inches	3-¾ inches
Average Anchor Embedment:	3-½ inches	3-½ inches
Minimum Anchor Embedment:	3 inches	3 inches
Available Alternate Anchor Locations:	Yes	Yes
Shipped Anchor Size:	8 inches	8 inches
Alternate Anchoring Methods:	Yes (see notes, below)	Yes (see notes, below)
Floor Levelness Requirement:	0.3125 (8mm) over 10 feet	0.3125 (8mm) over 10 feet

Table 6-2 Gantry and Table Mounting Requirements

- Note: If the Installers cannot set at least three anchors for the table, the installer must inform the customer that the minimum anchoring cannot be met, and a structural engineering contractor is required to determine the anchoring method and certify that their anchoring meets the stated GE minimum load requirement and torque specifications.
- Note: All other anchoring methods on floor types other than the concrete minimum must be determined at the customer's expense by a structural engineering contractor, and anchoring method must be certified to meet the stated GE minimum load requirement and torque specification.
- Note: If installing the GE LS scanner on a floor type other than a 4 inch concrete floor, all structural specifications in this document must be reviewed and met.

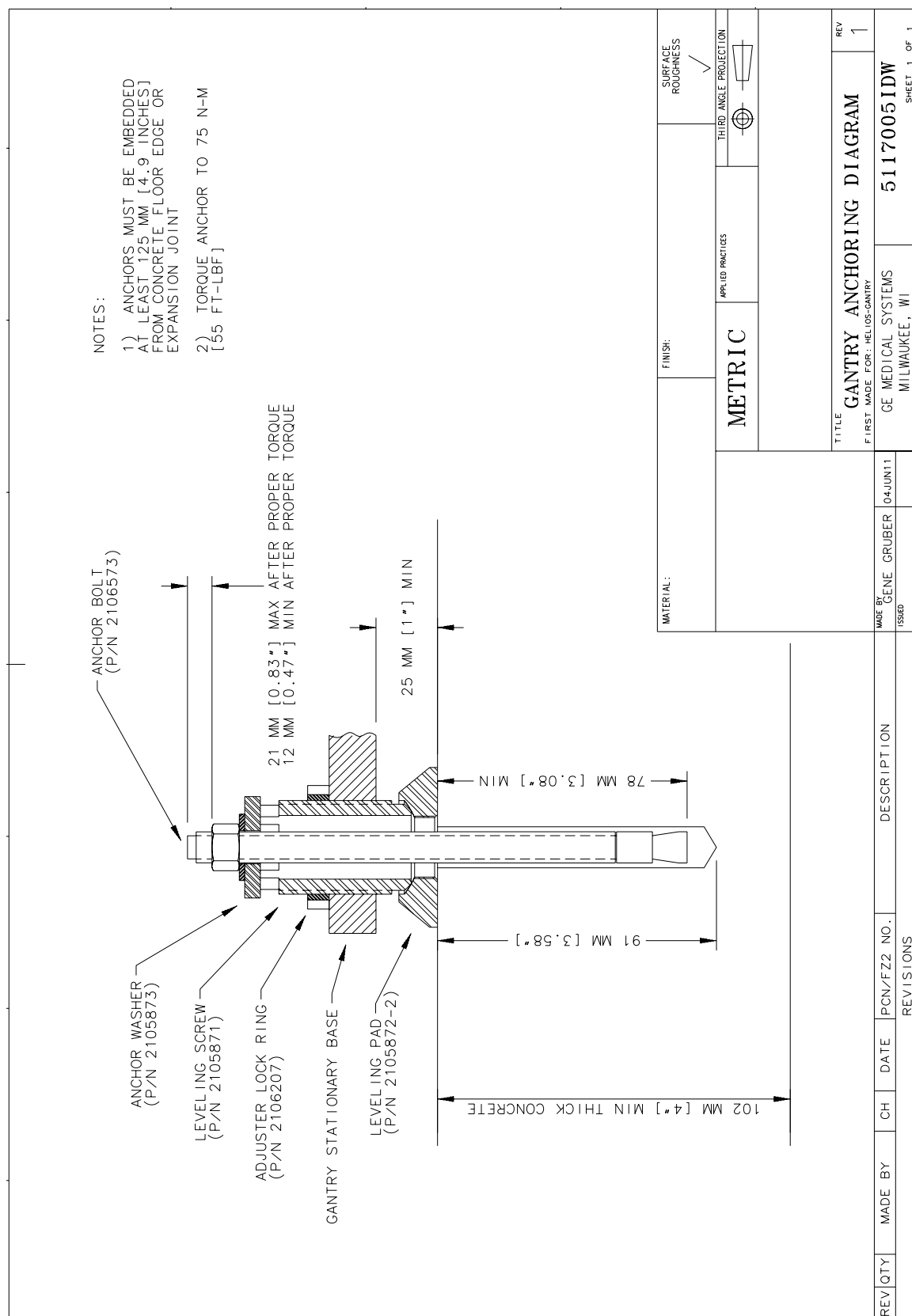


Figure 6-1 Typical Floor Anchor, Gantry and Table

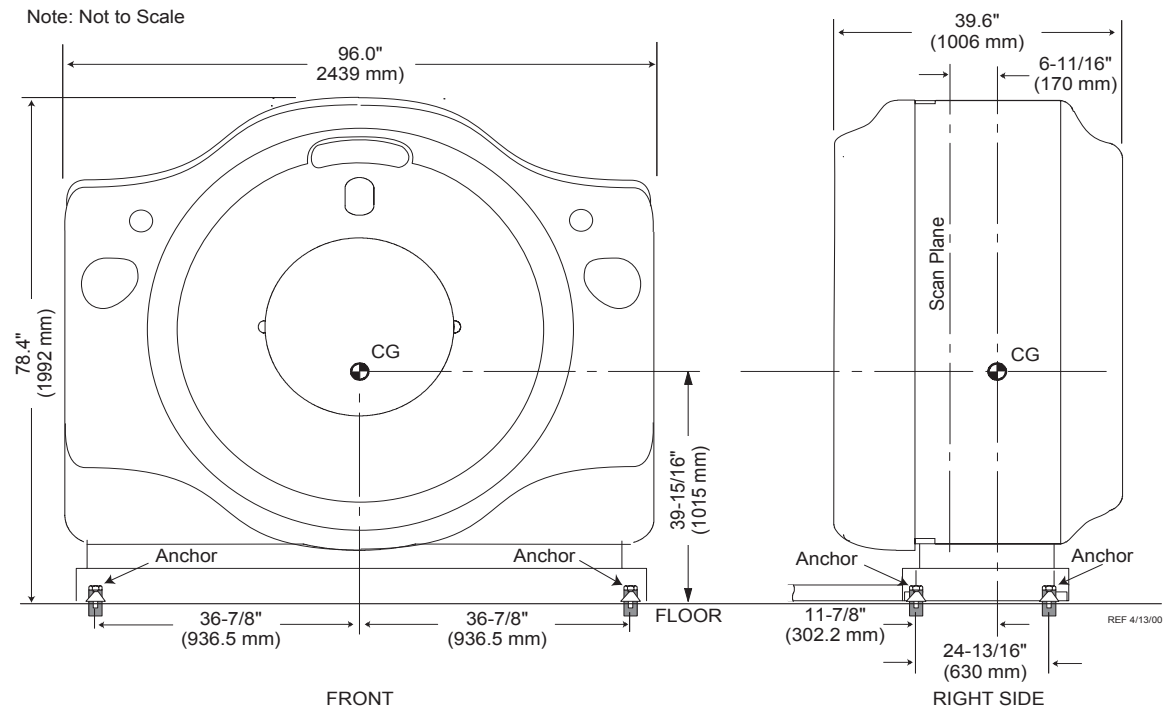


Figure 6-2 Gantry Center of Gravity

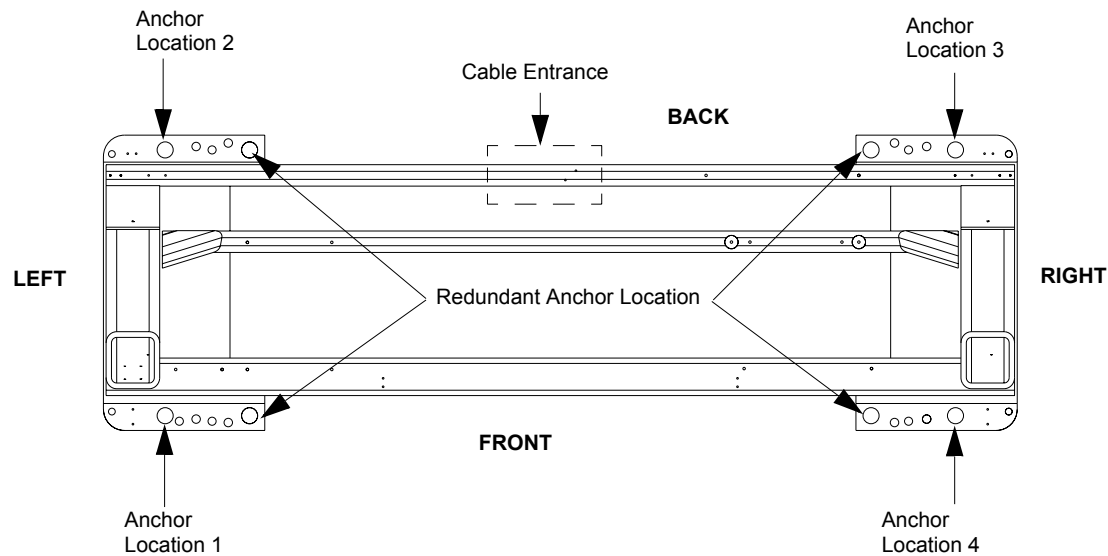


Figure 6-3 Gantry Anchor Locations

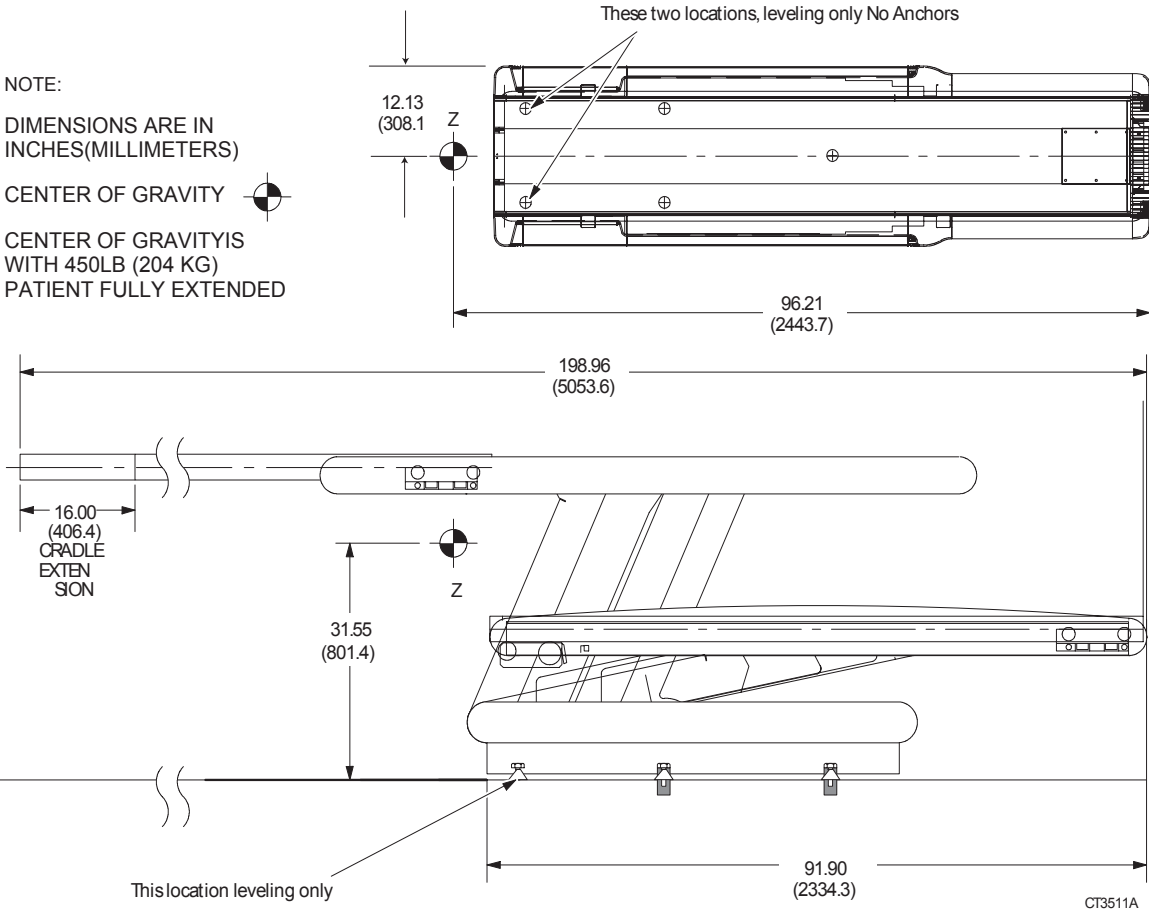


Figure 6-4 Patient Table (CT1)

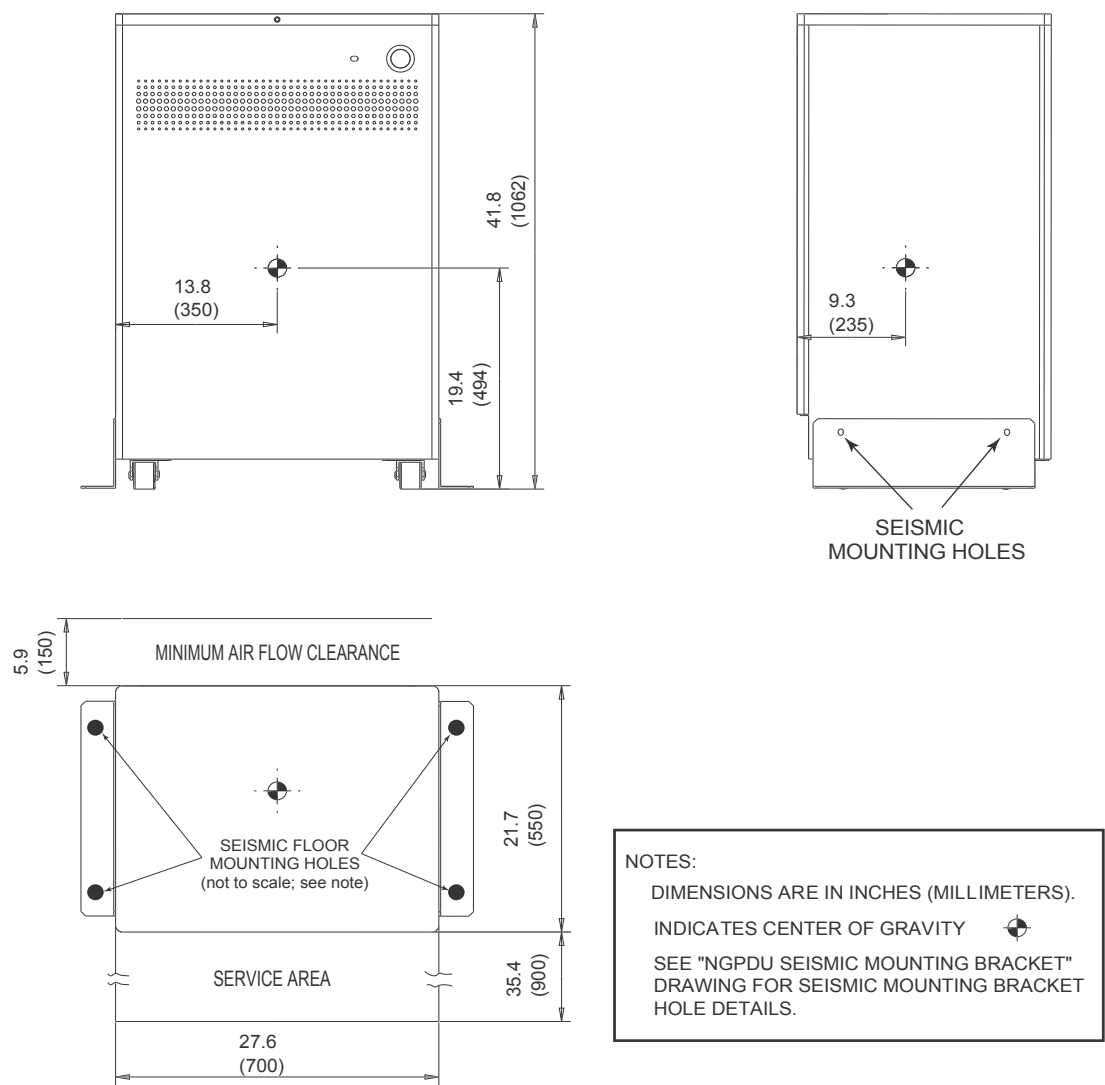


Figure 6-5 Power Distribution Unit (NGPDU)

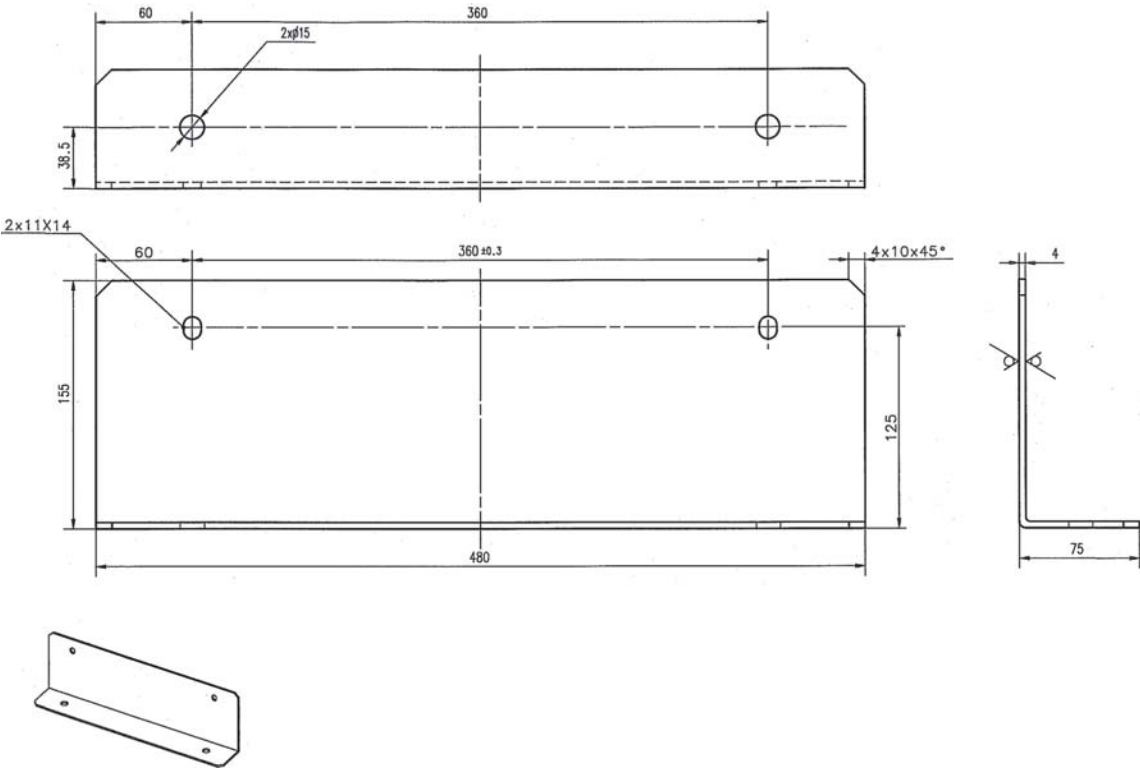


Figure 6-6 NGPDU Seismic Mounting Bracket (dimensions are in mm)

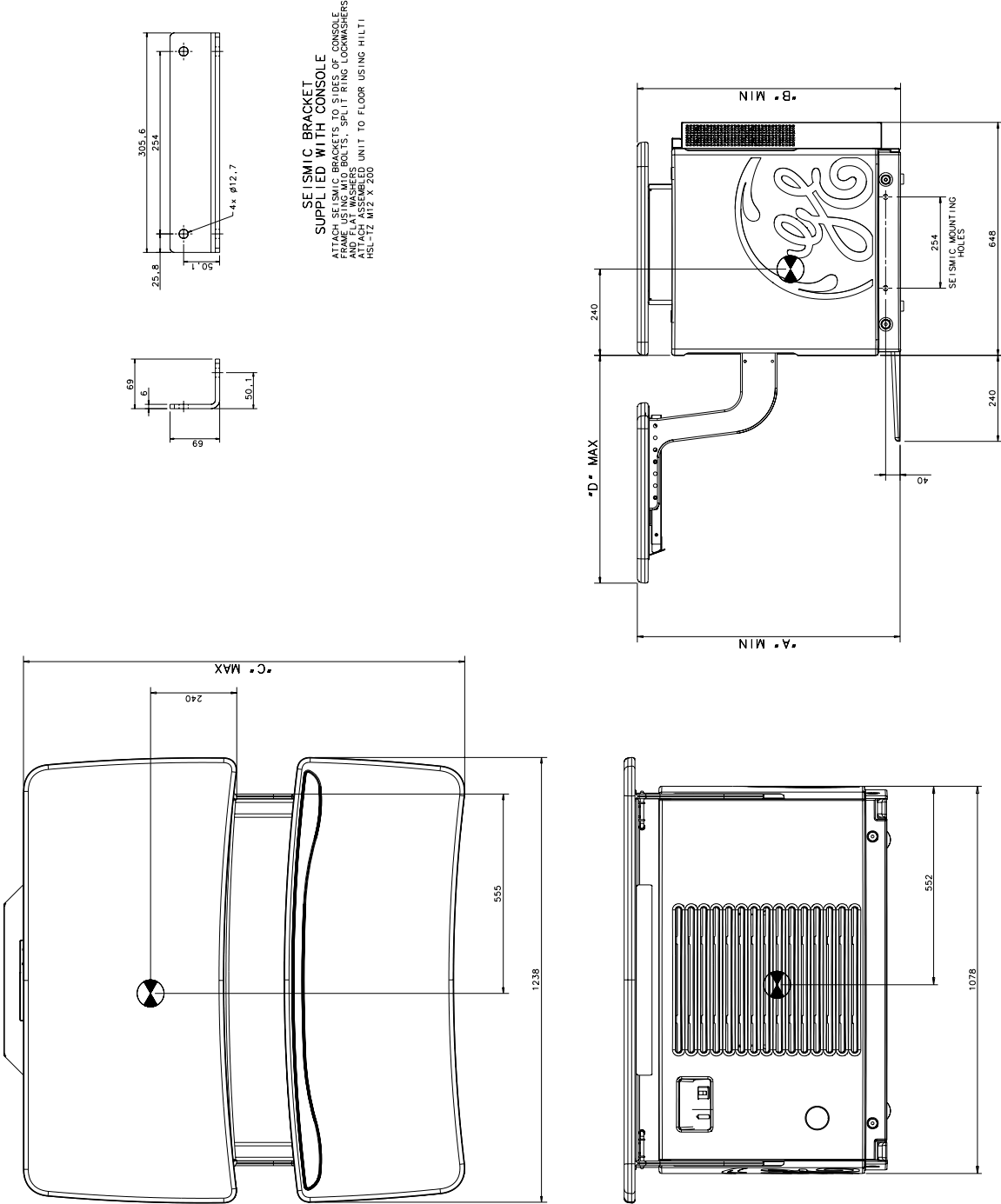


Figure 6-7 Operator's Console

Chapter 7

Delivery Data



WARNING SOME ASSEMBLIES ARE TOP-HEAVY. BE CAREFUL NOT TO TIP.

Section 1.0 Van Delivery

The CT system is packed for van shipment with minimum tear-down of components. It consists of approximately 20 shipping containers, which include dollies, skids and boxes without skids.

Section 2.0 Delivery/Shipping Requirements & Considerations

The LightSpeed RT System is not designed to tolerate excessive mishandling, including dropping, shock, vibration, tipping or hoisting.

The Gantry, Console, Table and PDU must NEVER be dropped. A drop from a height greater than 12.7 mm (0.5 in.) may induce structural damage to the frame or other major components. Damage resulting from a drop (e.g., bent frame, misalignment) may not be obvious until after system installation is complete.

To avoid dropping the Gantry, Dock to Dock shipment is recommended. Other methods are acceptable, provided that the system is not dropped or otherwise mishandled. For example, the system may be moved via flat-bed wrecker or by rolling it across SMOOTH sidewalks or other paved surfaces.

When moving Gantry off of a flat-bed wrecker, attach the straps to the lowest point possible on the dolly. Lower the Gantry at the slowest reasonable rate. See [Figure 7-1](#).

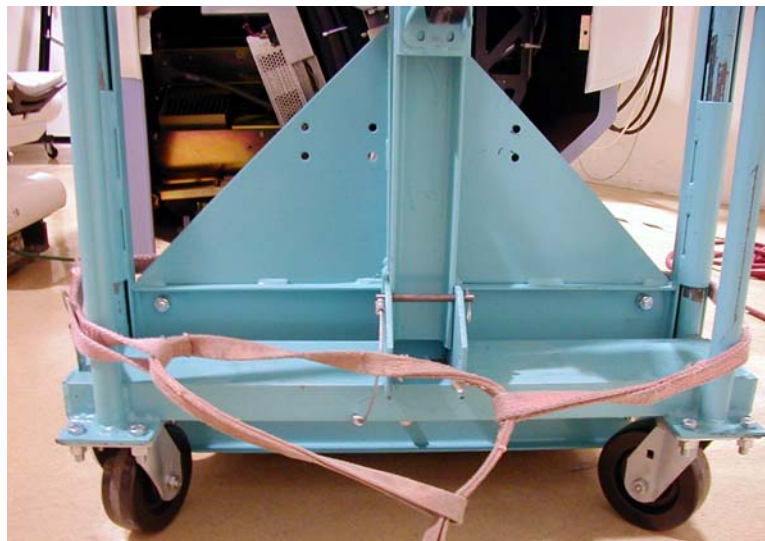


Figure 7-1 Gantry Strap Location

The LightSpeed RT System—including Gantry, Console, Table and PDU—is not designed to tolerate the excessive shock or vibration that may occur during unloading. For example, rolling the Console across a “washboard” style ramp may vibrate components to the extent of loosened or broken connections, etc. Damage resulting from shock or vibration (e.g., monitor, CD-ROM, hard-drive or octane failure) may not be evident until after system installation is complete.

All system components must remain upright at all times, and must not be tipped. Nor should the Gantry be hoisted. The Gantry should be moved by rolling it on its dollies only. Movement through hallways, doorways, elevators, etc., must be done without tipping or lifting the Gantry.

Protection for flooring along the move path (from dock to scan room) is advised.

**NOTICE**

Lifting a gantry requires engineering approval for each occurrence. Your Project Manager of Installation should contact CT Engineering for all special lifting requirements.

- **Do Not lift the gantry with a Forklift**
- **Unauthorized gantry lifting will cause Gantry Bearing Damage**

Engineer Contact Information:

Email: Yali.Zhang@geahk.ge.com

Phone: 86-10-58068355

Section 3.0

Site Environmental Considerations

3.1 Dust/Dirt Contamination

LightSpeed RT systems (consisting of: Console, PDU, Table and Gantry) are highly susceptible to airborne contaminants, especially concrete and drywall dust. Due to the possibility of contamination, these systems should NEVER be installed in a construction site. Any site with unfinished floors, walls or ceilings is considered a construction site, and is not suitable for system installation.

**NOTICE**

The act of installing a GE CT scanner in a construction (i.e., unfinished) site will likely result in the following adverse effects:

- **Increased installation time**
- **Decreased installation quality**
- **Increased scanner downtime, due to increased service calls**

3.2 Chemical Contamination

Wet film processors must never be installed in the same room as the scanner, due to the possibility of chemical contamination of LightSpeed RT CT scanner components. Such chemicals can contribute to increased equipment failures, increased system downtime, and decreased reliability. Film processor equipment installation must meet the manufacturer's requirements (e.g. ventilation specifications) and all applicable national and local codes. Also, consideration's should be given to the location of this equipment and chemical fumes relative to human contact as it relates to locating this equipment and chemicals in the control room.

Section 4.0

Storage Requirements



NOTICE Failure to adhere to Storage Requirements will likely result in equipment damage.

SHORT TERM STORAGE (LESS THAN 6 MONTHS)

If the CT system is to be stored before installation, store in a temperature and humidity controlled warehouse. Protect from weather, dirt and dust.

Meeting these requirements prevents rust and corrosion from forming on bearing surfaces due to condensation.

- Storage temperature should not exceed 4° to 27° C (40° to 80° F).
- Maintain relative humidity (non-condensing) between 20% and 60%.
- Maximum relative humidity rate of change is 5%/hr.
- The maximum temperature rate of change is 3° C/hr. (5° F/hr.)



NOTICE Between delivery is considered short-term storage.
Van storage must meet the same specifications as above.

LONG TERM STORAGE (6 MONTHS OR MORE)

If the CT system is to be stored before installation, store in a temperature and humidity controlled warehouse. Protect from weather, dirt and dust.

Meeting these requirements prevents rust and corrosion from forming on bearing surfaces due to condensation.

- Storage temperature should not exceed 10° to 32° C (50° to 90° F).
- Maintain relative humidity (non-condensing) between 20% and 70%.

Systems stored longer than six (6) months, but less than 12 months, should contact CT Engineering for installation start-up instructions. Storage longer than 12 months is not recommended.

Section 5.0

Extreme Temperature Transportation and Deliveries



NOTICE Failure to adhere to Extreme Temperature Transportation and Delivery requirements will likely result in equipment damage.

Extreme temperatures should be avoided during system transportation and delivery.

Extreme temperatures are defined as below -18° C (Zero degrees Fahrenheit) or above 49° C (120 degrees Fahrenheit), without humidity control.

Section 6.0

System Transportation

When transporting the CT system, ensure that the system is not exposed, for an extended period of time, to temperatures or humidity outside of the following specifications.

Temperature: -40° to +70° C (-40° to +158° F)

Humidity: 5% to 95%

**NOTICE**

Component Freezing occurs if CT system is exposed to temperatures below -18° C (0° F) for a period longer than two days.

Allow a minimum of 12 hours for the CT system to adjust to ambient room temperature, prior to installation.

Section 7.0

Gantry Considerations

The gantry is shipped with most covers installed. The assembly is mounted between two dollies. See [Figure 7-2](#). Two side rails are bolted to the dollies to stabilize dollies and protect gantry. Use dolly elevating casters to lift gantry off its base and roll it into position.



Figure 7-2 Gantry with Shipping Dollies and Side Rails

Door Openings. Clear door openings for moving equipment into building must be 1067mm X 2083mm (42 X 82 inches) minimum. 2439 mm (8 ft.) corridor width is helpful.

Elevator requirements. To move gantry from receiving location to scanning room, consider elevator capacity and size. By removing side rails and one dolly after gantry is placed in elevator, gantry width/length and elevator depth requirements are reduced. Contact a representative of elevator manufacturer if gantry weight exceeds elevator capacity.

CONFIGURATION	LENGTH	WIDTH	HEIGHT	WEIGHT
Dollies On, Side Rails On	3100 mm (122 in.)	1290 mm (51 in.)	2030 mm (79 in.)	2016 kg (4445 lb)
Dollies On, Side Rails Removed	3100 mm (122 in.)	1010 mm (39.4 in.)	2030 mm (79 in.)	1993 kg (4395 lb)
Dollies Off, Covers Off	2160 mm (85 in.)	860 mm (34 in.)	1870 mm (73 in.)	1655 kg (3648 lb)

Table 7-1 Size of Gantry & Dollies, with and without Side Rails

Minimum hallway and door size for gantry with covers and dollies attached, but side rails removed, is 1016 mm (40 in.).

For alternative lifting arrangements and instructions, contact GE Installation Support Services.

Dollies: Typically, dollies are used on the gantry, table, and console for domestic shipments. Once the gantry, table and console are installed at the site, return dollies to GE using the shipping document located in Box #1.

Dollies can be purchased for international shipments (B7850LD) to be used at the customer site. After the system has been removed from the crates, dollies shipped with international shipments

only are not to be shipped back to GE in Milwaukee, WI, USA, but to be retained to the local GE office or warehouse.

Zero clearance dollies are available (P/N 2128404). For information on zero clearance dollies, please go to: <http://www.umi-dollyshop.com/default1.htm>.

Section 8.0 Table Considerations

The table is shipped without side covers installed. Covers are shipped in four separate boxes. The table is mounted between two dollies.

Table shipping dimensions are 2565 mm (101 in.) long, 81 mm (32 in.) wide, and 119 mm (47 in.) high.

Section 9.0 Console Considerations

The console is shipped without the keyboard table installed. The keyboard table is shipped with the console.

- The console is shipped on a skid equipped with ramps for unloading.
- Do not remove the console from the shipping skid until it is in the CT equipment room.
- Console shipping dimensions (on the shipping skid) are 1016 mm (40 in.) deep, 1346 mm (53 in.) wide, and 1041 mm (41 in.) high.

Note: The dimensions of the operator console alone (as shipped) are 888 mm (34.95 in.) deep, 1238 mm (48.74 in.) wide, and 750 mm (29.5 in.) high.



NOTICE Do not lift the console by the monitor table top to remove from the shipping skid.



Figure 7-3 Console ready to be unloaded from shipping skid.

Chapter 8

Power Requirements

Section 1.0 Introduction

Note: Power requirements for LightSpeed RT differ from those for LightSpeed 4.X and previous systems. The power distribution unit (PDU) supplied with the LightSpeed RT system transforms and distributes power to all system components. The PDU is the only power entry point required to operate the system.

To minimize voltage regulation effects, power wiring between the facility main distribution panel and the PDU should be kept as short as possible.

When routing the power wiring all three phase wires and ground must be run in the same conduit or raceway duct. Power wires should be routed separately from system control and signal cables, using a separate conduit or trough in raceway duct.

Metallic conduit, floor duct or surface raceway may be used for running cables, depending upon local codes and practices. However, cable passageways should be large enough to install any cable with all other cables already installed. Use of non-metallic conduit is not recommended.

Section 2.0 System Input Power

2.1 Power Source Configuration

The LightSpeed RT CT Scanner is designed to operate on a three-phase, four-wire wye power source. A solidly grounded wye source is preferred. The neutral wire does not need to be run to the system, i.e., four-wire connection. If a neutral wire is run, then it should be terminated in the A1 box. A dedicated feeder from the nearest Main Distribution Panel (MDP) should supply power to the scanner. In accordance with the National Electric Code (U.S.) and similar applicable national and local codes, a protective disconnect device must be provided in the power line supplying the PDU. It must be located within 10 m (32 ft.) of the PDU, visible to PDU service personnel, and must have "lockout / tagout" provisions. This disconnect device is identified as "A1" in the interconnection schematic diagrams.

2.2 Rating

The LightSpeed RT system operates on three-phase power meeting the following specifications.

Voltage	380 to 480 VAC
Capacity	150 kVA
Frequency	50 or 60 Hz +/- 3 Hz

- Maximum power demand = 150 kVA @ 0.85 PF at a selected technique of 140 kV, 715 mA.
- Average (continuous) power demand at maximum duty cycle = 25 kVA.

The “A1” disconnect device referenced above must provide overcurrent protection for the system and facilitate multi-point remote “Emergency Off” control of system power. A disconnect utilizing undervoltage release control is preferred over shunt trip devices. The rating of the “A1” disconnect device depends on the nominal line voltage at the site. Refer to [Section 3.0: “Recommended Power Distribution System”](#), for minimum rating requirements and suggested disconnect devices.

2.3 Regulation

Total load regulation as measured at the PDU input terminals must not exceed 6%. The capacity of the facility transformer and size & length of feeder wires directly affect the load regulation presented to the system. Refer to [Section 3.0: “Recommended Power Distribution System”](#), for recommended single-unit installation specifics.

2.4 Phase Imbalance

The difference between the highest line-to-line voltage and lowest line-to-line voltage must not exceed 2% of the lowest line-to-line voltage.

2.5 Sags, Surges & Transients

Sags and surges of the power line must not exceed the absolute range limits shown in 103. Maximum transient voltages should be limited to 1500V peak.

2.6 Grounding

Metal conduit, raceway or the armor of armored cable used to power the system should be bonded to the PDU cabinet. However, in addition to such mechanical grounding, a dedicated 1/0 (55 mm²) or larger insulated copper ground wire must be run with the phase wires from the main distribution panel to the PDU.

Note: The shield or armor of armored cable is not sufficient for this purpose.

The ground wire should be bonded to intermediate distribution panels through which it passes in accordance with local codes. The resistance between the PDU ground and the facility earth ground must not exceed 0.5 ohm. In addition, the total resistance between the PDU ground and earth must not exceed 2 ohms.

Section 3.0

Recommended Power Distribution System

A dedicated feeder run from the facility main isolation transformer is recommended to power the LightSpeed RT CT scanner. If the scanner must be powered from an existing distribution transformer and secondary feeder, such as the equipment distribution panel of an X-ray department, installation with other X-Ray equipment that use rapid film changers should be avoided. These changers use a large number of high powered, closely spaced exposures, which may coincide with the CT scan and produce image artifacts.

If a dedicated distribution transformer is provided for the scanner, the minimum recommended transformer size is 225 kVA, rated 2.4% regulation at unity power factor. For this configuration, the minimum recommended feeder size and overcurrent protection device based on line voltage is shown in 103 Minimum Feeder Wire Size.

In all cases, qualified personnel must verify that the transformer and feeder, at point of take-off, plus the run to the LightSpeed RT CT scanner meet all the requirements stated in this document.

SYSTEM CHARACTERISTICS:

- Maximum power demand = 150kVA @ 0.85 PF: at a Selected Technique of 140 kV, 715 mA.
- Continuous (average) power demand at maximum duty cycle = 25kVA.
- Maximum allowable total source regulation is 6%.
- Minimum recommended transformer size: 225 kVA, with 2.4% rated regulation at unity power factor. Resultant maximum allowable feeder regulation is 3.4%.

The nominal line voltage must fall within one of the ranges listed below

Nominal Line Voltage	380	400	420	440	460	480
Hi-Line Limit, +10%	418	440	462	484	506	528
Lo-Line Limit, -10%	342	360	378	396	414	432
Continuous Line Current	38	36	34	33	31	30
Momentary Line Current	228	217	206	197	188	180
Maximum Line Current	253	241	229	219	209	200
Minimum Recommended Circuit Protection Rating	150	150	150	125	125	125

Table 8-1 Nominal Line Voltage

OPTION	CAT NUM
225 kVA Transformer	E4500AW
Isolation Transformer	E4500BC

Table 8-2 Purchasable Options

FEEDER LENGTH (MDA TO A1) FEET (METERS)	MINIMUM FEEDER WIRE SIZE, AWG OR MCM (SQ. MM)/ VAC					
	380 VAC	400 VAC	420 VAC	440VAC	460VAC	480VAC
50 (15)	1/0 (55)	1/0 (55)	1/0 (55)	1 (45)	1 (45)	1 (45)
100 (30)	1/0 (55)	1/0 (55)	1/0 (55)	1 (45)	1 (45)	1 (45)
150 (46)	1/0 (55)	1/0 (55)	1/0 (55)	1 (45)	1 (45)	1 (45)
200 (61)	1/0 (55)	1/0 (55)	1/0 (55)	1 (45)	1 (45)	1 (45)
250 (76)	2/0 (70)	2/0 (70)	1/0 (55)	1/0 (55)	1 (45)	1 (45)
300 (91)	3/0 (85)	3/0 (85)	2/0 (70)	2/0 (70)	1/0 (55)	1/0 (55)
350 (107)	4/0 (100)	3/0 (85)	3/0 (85)	2/0 (70)	2/0 (70)	1/0 (55)
400 (122)	250 (125)	4/0 (100)	3/0 (85)	3/0 (85)	3/0 (85)	2/0 (70)

Table 8-3 Minimum Feeder Wire Size

SUB-FEEDER LENGTH (A1 TO PM) FEET (METERS)	MINIMUM SUB-FEEDER WIRE, AWG OR MCM (SQ. MM)					
	380 VAC	400 VAC	420 VAC	440VAC	460VAC	480VAC
32 (9.7536)	1/0 (55)	1/0 (55)	1/0 (55)	1 (45)	1 (45)	1 (45)

Table 8-4 Minimum Sub-Feeder Wire Size

- 1.) 103, 103, and 103 above are based on the use of copper wire, rated 75 C and run in steel conduit. Ampacity is determined in accordance with the National Electrical Code (NFPA 70), Table 310-16 (2002)
- 2.) The minimum feeder size is determined by the ampacity of the circuit protection device listed above, except where a larger size is necessary to meet total source regulation limits.
- 3.) A 1/0 (55 sq. mm) ground wire is recommended in all cases.

Section 4.0 Ground System

The LightSpeed RT CT Scanner has been designed to use an equal potential grounding system. The required ground system is shown in [Figure 8-1: "System Ground Map"](#). There are three primary grounding points:

- A system power ground point located in the PDU.
- A reference ground point located between gantry and table base.
- A patient ground point located at the front of the table base.

All exposed metal surfaces in the patient vicinity are grounded to the reference ground point.

For additional information, refer to Electrical Safety Equipment, Direction 46-014505.

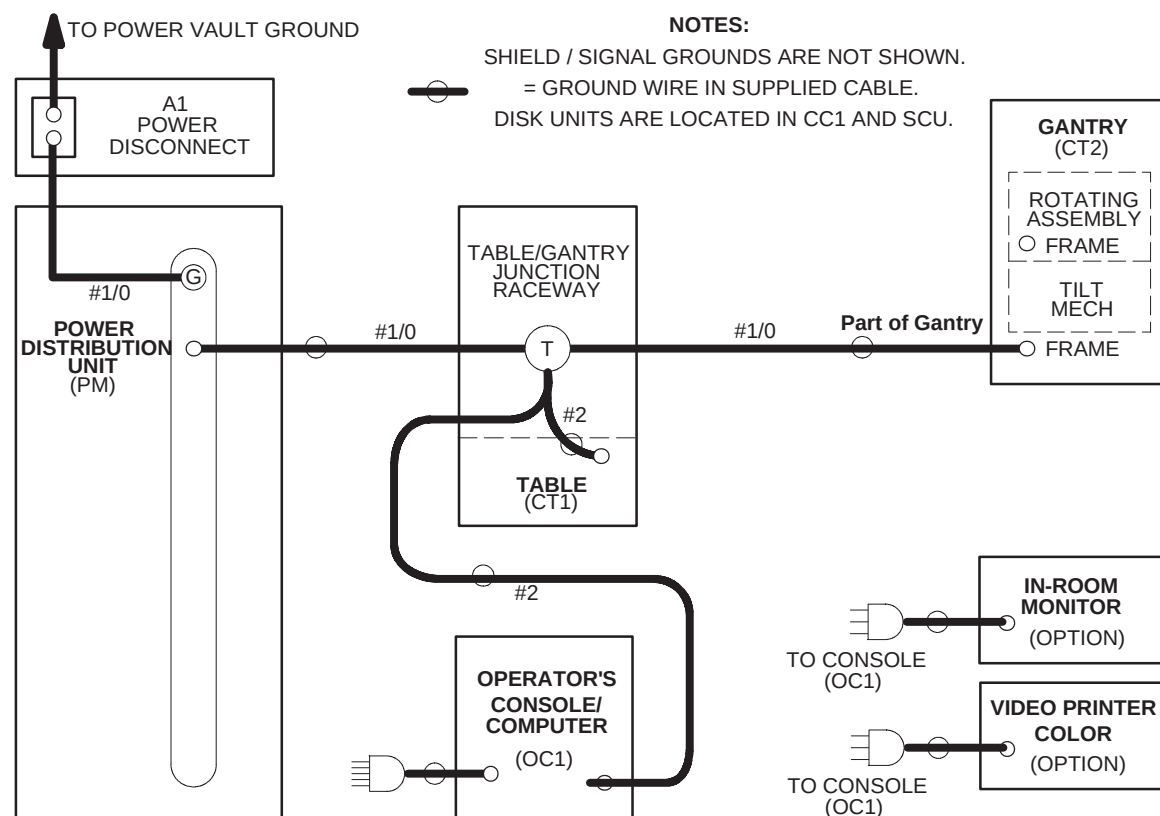


Figure 8-1 System Ground Map

This Page Intentionally Blank.

Chapter 9

Interconnection Data

Section 1.0 Introduction

Figure 9-2 shows interconnection runs for a 50/60 Hz system.

Table 9-1 shows component designators for supplied equipment and options and wall power outlets.

Table 9-5 lists customer-installed wiring and supplied cables. Actual length of each run is less than the length of supplied cables to allow for routing inside equipment. Cable diameters and sizes of connectors are provided to aid in sizing conduit and access plates.

Table 9-2 and Table 9-3 list details for connection to LightSpeed Series equipment, using standard (short) length and non-standard (long) length cables, respectively. Details are listed for the following types of runs as appropriate:

- Flush-floor duct
- Computer floor
- Through-wall bushing
- Junction box
- Surface floor duct
- Through-floor duct
- Wall duct
- Conduit

Need for additional junction boxes is minimized by use of either a cable raceway system or a raised computer floor. LightSpeed Series systems use prefabricated cables with large plugs. Therefore, conduit or pipe is not recommended for cable runs.

Note: Use dry cleaning for electro components.

Section 2.0 Component Designators

DESIGNATOR	APPLIES TO	SOURCE
A1	Primary power disconnect	Contractor supplied
CT1	Patient table	System
CT2	Gantry	System
ITL	InSite telephone lines	Contractor supplied
OC1	Operator's console/computer	System
PDU	Power distribution unit	System
SEO	System emergency off	Contractor supplied
SM	Slave monitor	Option
WL	"X-ray on" warning light	Contractor supplied
DS	Door Interlock Switch	Contractor supplied
BBNC	Broad-Band Network Connection	Contractor supplied

Table 9-1 Component Designators

Section 3.0 Interconnect Runs, Wiring and Cables

3.1 GE Healthcare Supplied (Standard Length)

RUN #	LENGTH, ACTUAL (USABLE)		PART #	DESCRIPTION	UL CABLE INFORMATION								PULL SIZE MM (INCHES)
	ft	m			UL Style	Flam. Rating	Voltage Rating	Voltage Actual	Temp. Rating (C)	Dia. mm (inch)	# of Cond	Size AWG	
050	28 (20)	8.5 (6.1)	2343529-2	HVDC, PDU to Gantry	2587	FT4	600	+ & - 350VDC	90	19 (.751)	3	(2) 4 (1) 8	22 (.87) Dia
051	28 (20)	8.5 (6.1)	2343530-2	HVAC, PDU to Gantry	2587	FT4	600	440Y/254	90	15.3 (.604)	4	14	11.2 (.44) Dia
052	28 (20)	8.5 (6.1)	2343528-2	LVAC, PDU to Gantry	2587	FT4	600	208Y/120	90	13.8 (.542)	5	8	56.4 (2.22) Dia
053	65 (60)	19.8 (18.3)	2343531-2	LVAC, PDU to Console	2587	FT4	600	120VAC	90	12.3 (.483)	3	10	56.4 (2.22) Dia
054			2343596	LVAC, Gantry to Table	1015		600	120VAC			3	14	
055	43 (35)	13.2 (10.67)	2371450-2	Ground, PDU to Raceway	1284	VW-1 (FT-1)	600	0	105	15.5 (.608)	1	1/0	15.8 (.62) Dia
056	71 (60)	21.7 (18.3)	2371450-4	Ground, Raceway to Console	1283	VW-1 (FT-1)	600	0	105	11.9 (.467)	1	2	12.2 (.48) Dia
100	32.5 (20)	9.9 (6.1)	5120646-2	Signal, Gantry MSUB to PDU		FT-4	300	<30VDC	80	11.2 (.440)	25	22	17 x 58 (.68 x 2.30) 19 x 51 (.75 x 2.01)
101	71 (60)	21.7 (18.3)	5120645-2	Signal, Gantry MSUB to OC		FT-4	300	<30VDC	80	11.2 (.440)	25	22	17 x 58 (.68 x 2.30) 19 x 51 (.75 x 2.01)
102	68 (60)	20.7 (18.3)	2352714-3	Signal (LAN), Gantry to OC			1900	<30VDC		5.9 (.234)	8	24	15 (.59) Dia
103	68 (60)	20.7 (18.3)	2117848-7	Fiber Optic, Gantry to OC			N/A	N/A			1	N/A	10 (.39) Dia
104			2333151	Signal, Gantry to Table		FT-4	300		80		25	22	

Table 9-2 GE Healthcare Supplied Cables (Standard Run) - UL Information

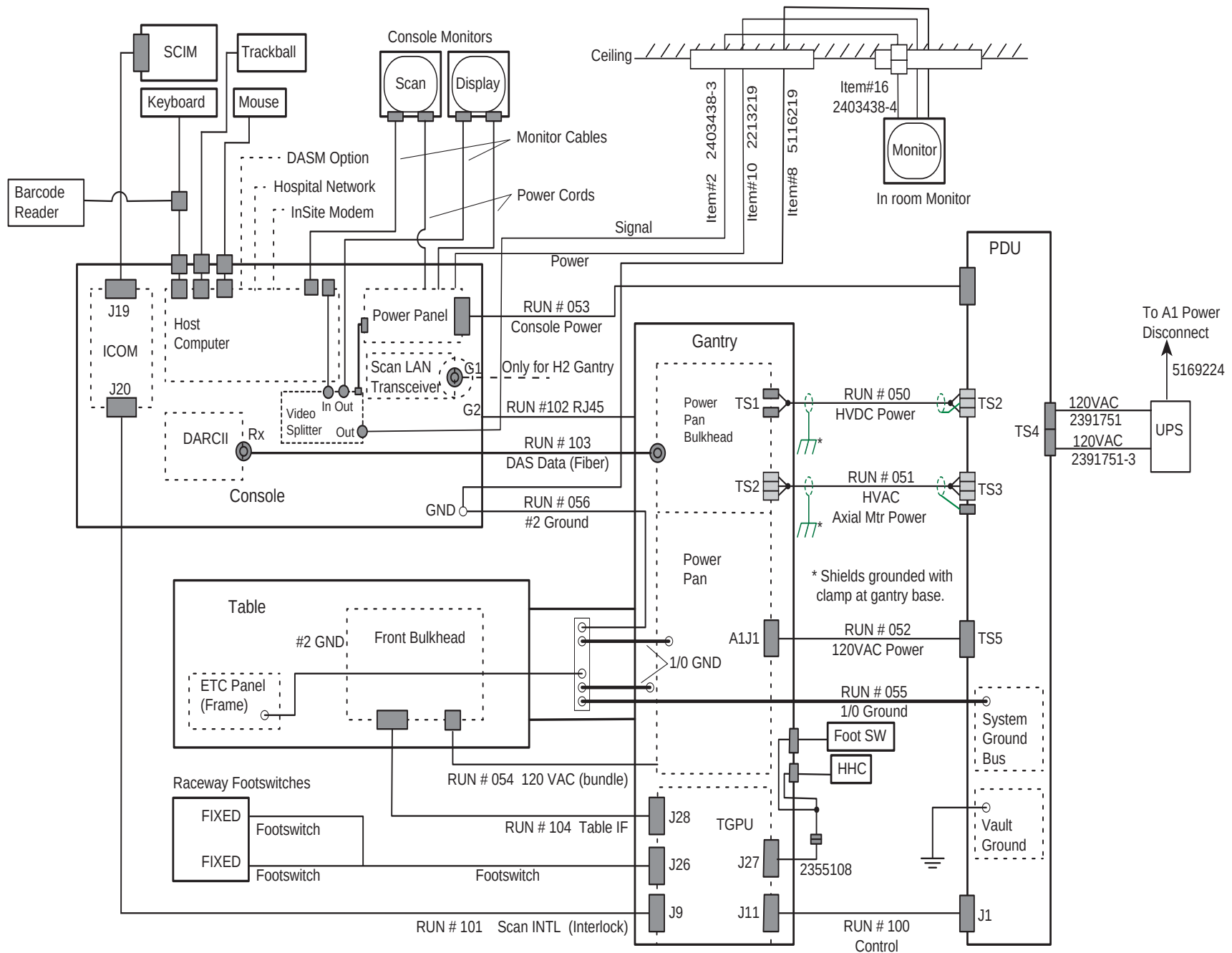


Figure 9-1 System Interconnect Diagram

3.2 GE Healthcare Supplied (Optional, Long Run)

RUN #	LENGTH, ACTUAL (USABLE)		PART #	DESCRIPTION	UL CABLE INFORMATION								PULL SIZE MM (INCHES)
	ft	m			UL Style	Flam. Rating	Voltage Rating	Voltage Actual	Temp. Rating (C)	Dia. mm (inch)	# of Cond	Size AWG	
050	63 (55)	19.3 (16.76)	2343529	HVDC, PDU to Gantry	2587	FT4	600	+ & - 350VDC	90	19 (.751)	3	(2) 4 (1) 8	22 (.87) Dia
051	62.5 (55)	19 (16.76)	2343530	HVAC, PDU to Gantry	2587	FT4	600	440Y/254	90	15.3 (.604)	4	14	11.2 (.44) Dia
052	60 (55)	18.5 (16.76)	2343528	LVAC, PDU to Gantry	2587	FT4	600	208Y/120	90	13.8 (.542)	5	8	56.4 (2.22) Dia
053	80 (75)	24.5 (22.86)	2343531	LVAC, PDU to Console	2587	FT4	600	120VAC	90	12.3 (.483)	3	10	56.4 (2.22) Dia
055	63 (55)	19.3 (16.76)	2371450	Ground, PDU to Raceway	1284	VW-1 (FT-1)	600	0	105	15.5 (.608)	1	1/0	15.8 (.62) Dia
056	83 (75)	25.5 (22.86)	2371450-3	Ground, Raceway to Console	1283	VW-1 (FT-1)	600	0	105	11.9 (.467)	1	2	12.2 (.48) Dia
100	63 (55)	19.3 (16.76)	5120646	Signal, Gantry MSUB to PDU		FT-4	300	<30VDC	80	11.2 (.440)	25	22	17 x 58 (.68 x 2.30) 19 x 51 (.75 x 2.01)
101	83 (75)	25.5 (22.86)	5120645	Signal, Gantry MSUB to OC		FT-4	300	<30VDC	80	11.2 (.440)	25	22	17 x 58 (.68 x 2.30) 19 x 51 (.75 x 2.01)
102	80 (75)	24.3 (22.86)	2352714-2	Signal (LAN), Gantry to OC			1900	<30VDC		5.9 (.234)	8	24	15 (.59) Dia
103	80 (75)	24.3 (22.86)	2117848-2	Fiber Optic, Gantry to OC			N/A	N/A			1	N/A	10 (.39) Dia

Table 9-3 GE Healthcare Supplied Cables (Optional, Long Run) - UL Information

3.3 GE Healthcare Supplied (Cables of Options)


OPTION	LENGTH, ACTUAL (USABLE)		PART #	DESCRIPTION	UL CABLE INFORMATION								PULL SIZE MM (INCHES)
	ft	m			UL Style	Flam. Rating	Voltage Rating	Voltage Actual	Temp. Rating (C)	Dia. mm (inch)	# of Cond	Size AWG	
Fluoro	75	22.9	2403438-3	5 BNC MALE TO HD 15 MALE 75 FEET	1015	FT4		1Vp-p	75	9.1 (0.358)	5	26	
	70	21.2	2213219	POWER CABLE FOR LCD-CONSOLE TO LCD		FT1	120	120VAC	105	9.3 (0.366)	3	14	
	71	21.5	5116219	Grounding Cable For LCD Console To LCD		VW-1	600	0V	105		1	8	
	15	4.6	2403438-4	HD 15 FEMALE TO HD 15 MALE 15 FEET				1Vp-p	60	8.0 (0.315)	5	26	
	1.3	0.4	2355108	JUMPER CABLE FOR ADAPTING 2286150 TO WORK WITH H-POWER MSUB							8	22	
UPS	15	4.6	2391751	POWER CABLE, NGPDU TO UPS	2587	FT4	600	208VAC	90	5.8 (0.228)	5	8	
	15	4.6	2391751-3	POWER CABLE, UPS DISCONNECT PANEL TO NGPDU	2587	FT4	600	208VAC	90	5.8 (0.228)	4	8	
	45	13.6	5169224	UPS CONTROL CABLE	2587	FT4	600	120VAC	90	10.3 (0.406)	5	18	

Table 9-4 GE Healthcare Supplied Cables for Options - UL Information

3.4 Contractor (Customer) Supplied

CUSTOMER INSTALLED WIRING		DESCRIPTION	CABLES SUPPLIED			PLUG PULLING DIMENSIONS		WIRE & CABLE PIGTAILS FT. (M.)	
QTY	SIZE AWG (MM ²)		PART NO	LENGTH FT. (M.)	DIA. IN (MM)	FROM	TO	FROM	TO
RUN NO. 1 FROM PRIMARY POWER SOURCE TO FACILITY DISCONNECT (POWER SOURCE - A1)									
Maximum Run Length *									
3	*	POWER						3 (1)	3(1)
1	1/0 (50)	GROUND						3 (1)	3 (1)
RUN NO. 2 FROM FACILITY DISCONNECT TO POWER MODULE (A1 - PM) MAXIMUM RUN LENGTH *									
3	*	POWER						3 (1)	3(1)
1	1/0 (50)	GROUND						3 (1)	3 (1)
RUN NO. 3 FROM FACILITY DISCONNECT TO SYSTEM EMERGENCY OFF (A1 - SEO)									
2	14 (2)	POWER						6 (2)	6 (2)
1	14 (2)	GROUND						6 (2)	6 (2)
RUN NO. 4 POWER MODULE TO WARNING LIGHT CONTROL (PM - WL)									
2	14 (2)	WARNING LIGHT 24 VOLT CONTROL A3J2-1,2,3,4							
RUN NO. 5 POWER MODULE TO SCAN ROOM DOOR INTERLOCK (PM - DOOR SWITCH)									
2	14 (2)	SCAN ROOM DOOR INTER LOCK A3J6-1,2							
*	REFER TO Table 8-3 on page 103 FOR AWG (MM2) WIRE SIZES.								

Table 9-5 Runs 1, 2, 3, 4 and 5 Connections

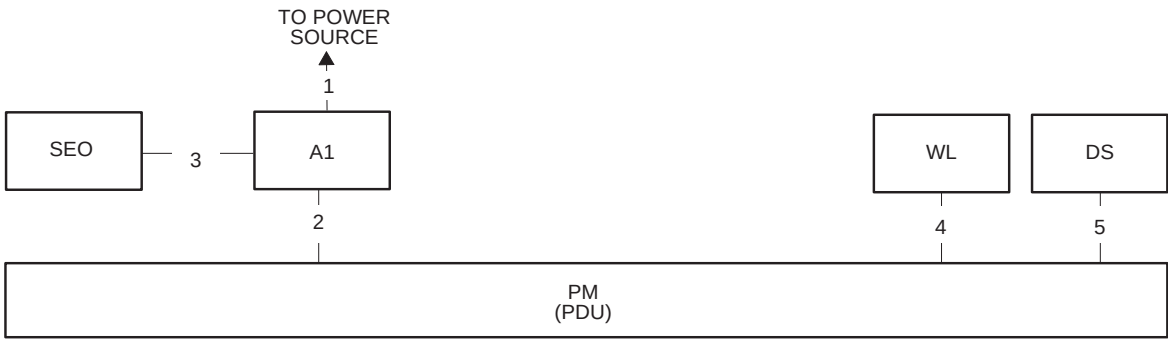

NOTES:

- 1) Used for remote diagnostics - Option
- 2) Refer to the appropriate Pre-installation / Installation documents for the Laser Camera
- 3) Category 5 cable. Use one of the following patch cords:

CAT Num	GE Part Num	Length
K9000WB	2215028-10	20 m
K9000KP	2215028-5	10 m
K9000JR	2215028-4	5 m
K9000WA	2215028-9	3 m

- 4) In order to avoid any violation of each National Regulation (NEC in USA, CCC in China, etc.), use of the compliant cable/wire is recommended. For the China market, the China end-user shall purchase the power supply cable that has the CCC mark.

Only one phone connection is required for the system.

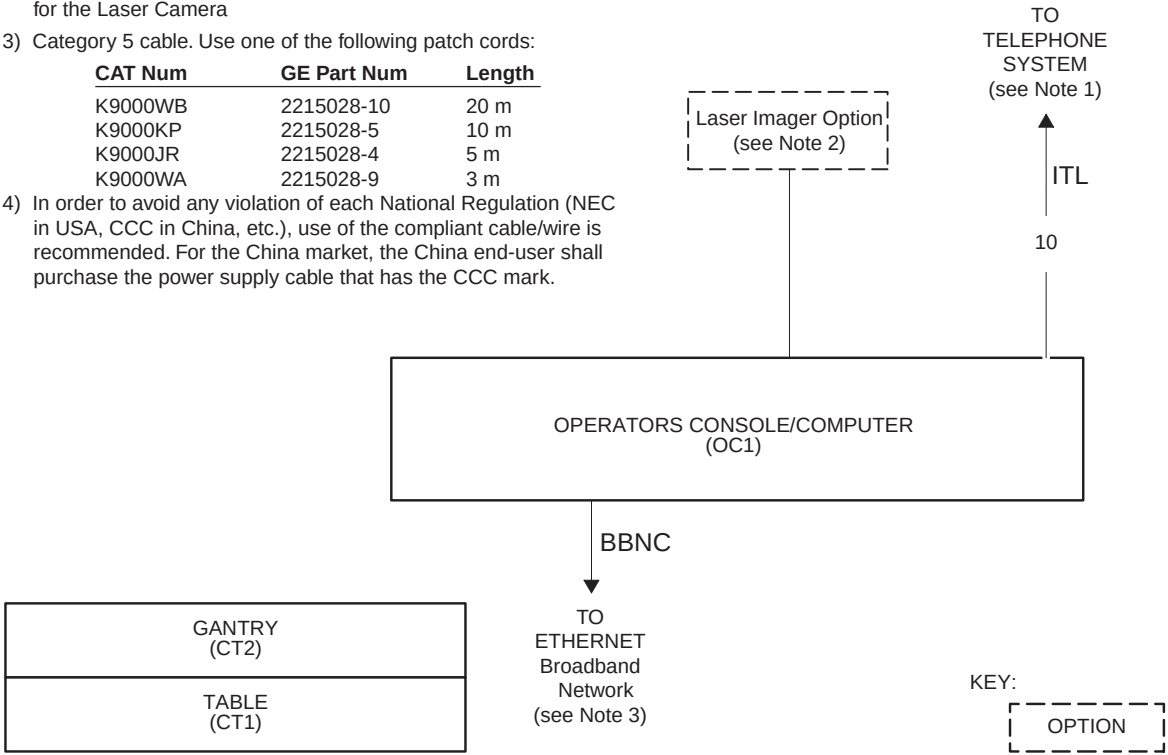


Figure 9-2 Interconnection Runs

Section 4.0

Contractor Supplied Components

REFERENCE	ASSOCIATED EQUIPMENT	MATERIAL/LABOR SUPPLIED BY CUSTOMER CONTRACTOR	USA VENDOR / CAT NO. GE CATALOG
A1 380 - 480V 50/60 Hz	Circuit Breaker with Magnetic Contactor	3 Pole, 380V - 480V, Combination breaker with magnetic contactor. Includes control transformer, optional UPS interface, On/Off controls and auto-restart feature	Recommend*: - E4502AE (125A) - E4502AF (150A) Optional remote operator control available from GE Supply, Cat # GESCTR0CS1
BBNC (required)	Broad-Band Network Connection	Broad-Band network connection wall jack, located within 1m (39inches) of console location, for internal hospital networking and InSite Broad-Band connectivity. Cabling to conform to facility's IT standards.	
ITL (optional)	In-suite Telephone Lines	Supply 2 voice-grade telephone lines. One line must be a direct number from outside the facility – do not route this line through a telephone switchboard. Telephone line operating charges are paid by customer.	
	System Components	Reference the system installation drawings supplied by Installation Support Services within your geographic area.	

*Refer to Table 9-8 on page 116.

Table 9-6 Contractor-Supplied Components

Section 5.0

Fuse

MARK	NUMBER	QTY	FRU	DESCRIPTION/NAME
1	2351492	2.000	Yes	200A FUSE
2	2364059	2.000	Yes	GLASS FUSE
3	46-170021P50	2.000	Yes	FUSE 12 AMPS 250 VOLTS
4	46-170021P106	2.000	Yes	FUSE 8 AMPS 250 VOLTS
5	2336517-2	2.000	Yes	FUSE 25 AMPS 700 VOLTS 2.
6	46-170021P52	3.000	Yes	FUSE 3 AMPS 250 VOLTS
7	46-170021P10	2.000	Yes	006.000A 0250V 3AG
8	46-170021P101	1.000	Yes	FUSE 20 AMPS 700 VOLTS
9	46-170021P11	2.000	Yes	7.0A, 250V SLO-BLO FUSE.
10	46-170021P31	2.000	Yes	1/2A, 250V SLO-BLO FUSE.
11	2106993-5	2.000	Yes	FUSE 20 AMPS 500 VOLTS .3
12	46-327160P1	2.000	Yes	12.0A, 125VAC, DUAL ELEME
13	2379651	1.000	Yes	FUSE - 700 V, 200 A, FAST

Table 9-7 H-Power Fuse Kit (2218570-2 BOM, rev 2)

Section 6.0

UPS Interconnect

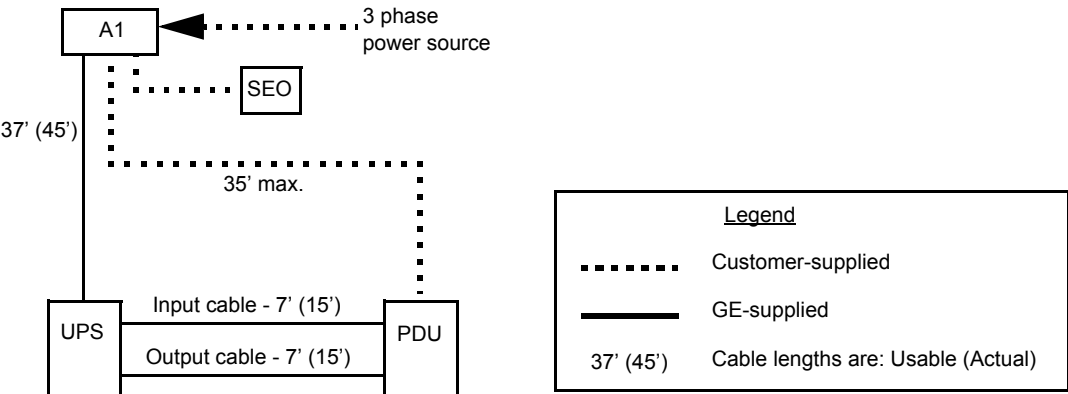


Figure 9-3 Typical UPS Interconnect

The UPS Kit catalog number is B7999ZA. See table below, for detail.

NOTICE **UPS Kit B7999ZA REQUIRES** installation of one of the A1 Panels listed below.

PDU Type & Model #	Max. Mom. kVA Rating	Recommended Main Disconnect (A1) Cat #		Optional Partial UPS Kit Cat #
		Europe & Asia (380-400V or 420V)	North America (440V or 460-480V)	
NGPDU 2326492	150kVA	E4502AF (150A) (incl. Auto Restart & Integrated UPS Control)	E4502AE (125A) (incl. Auto Restart & Integrated UPS Control)	B7999ZA alt. E4502KY (Includes PW 9155- 10GE10kVa, 2ph. UPS,& hardware kit) REQUIRES one of the A1 Panels shown at left, or equal.

Table 9-8 LightSpeed RT Partial UPS Back-up Options

Conduit is required between:

- A1 and UPS
- UPS and PDU
- PDU and A1

Section 7.0


Typical Customer Supplied Wiring

7.1 Primary Power Disconnect

Installing this system requires a Lockout/Tagout-compatible disconnect. If a UPS is required, a GE Disconnect is strongly recommended for safe operation. The GE disconnect and UPS are designed to work together.



Figure 9-4 Primary Power Disconnect (A1) - Typical GE A1 shown

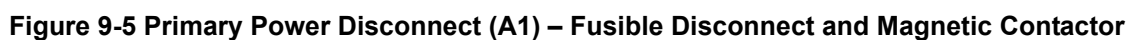

7.2 Scan Room Warning Light & Door Interlock

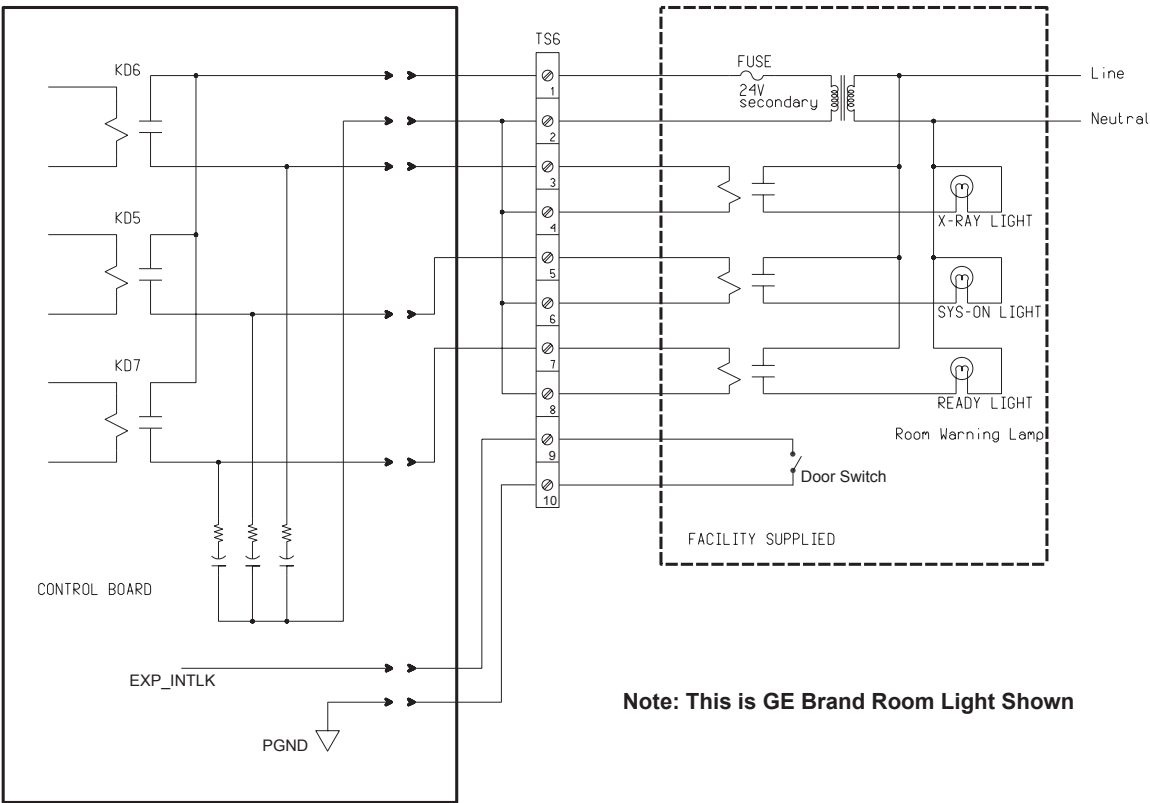


Figure 9-6 TS6 X-Ray Warning Light Connections

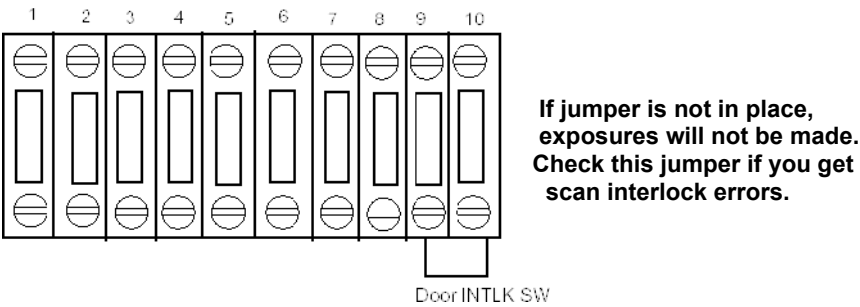


Figure 9-7 TS6 Room Door Interlock Connections - Without a Door Interlock

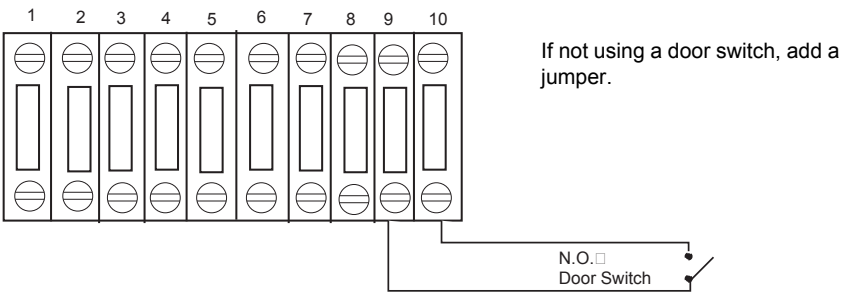


Figure 9-8 TS6 Room Door Interlock Connections - With a Door Interlock

© 2003-2006, General Electric Company.
GE Medical Systems, a General Electric Company, going to market as GE Healthcare.
3000 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com

