



GE Medical Systems

Technical Publication

**Direction 2249157FMI
Revision 0**

GE Medical Systems FMI 25258 2-Way Video Splitter Replacement

Restricted To General Electric Medical Systems (GEMS)

This Material is the property of GEMS and contains information which is restricted to use by GEMS Manufacturing, Engineering and Service personnel. Any reproduction in whole or in part is strictly prohibited. Neither this material nor any information contained herein shall be disclosed to anyone outside of GEMS Manufacturing, Engineering and Service personnel.

Copyright ©1999 by General Electric Company



Restricted to General Electric Medical Systems (GEMS) Personnel Only.

Permission to use this Restricted Service Software and related documentation (hereinafter called the "Material") by persons other than GE Medical Systems employees is prohibited. This Material is the property of GEMS and contains information which is restricted to use by GEMS Manufacturing, Engineering and Service personnel. Any reproduction in whole or in part is strictly prohibited. Neither this material nor any information contained herein shall be disclosed to anyone outside of GEMS Manufacturing, Engineering and Service personnel.

Revision History

Rev	Date	Reason for Change
Rev. 0	Date: 06/28/99	Initial Release

List of effective pages

Page(s)	Rev.
Title Page	0
Restricted/History/LOEP	0
Pages 1 to 5	0
Rear	-

2-Way Video Splitter Replacement

EFFECTIVITY:

All 2-way splitters from catalog # B7530RB & RC between 6/23/97 to 3/10/99.

PREREQUISITE:

None

PURPOSE:

To replace the current two-way video splitter (part # 2151622) for remote monitor options.

MANPOWER REQUIREMENTS:

One FE needed for no more than .5 hours labor and no more than 2 hours travel.

SPECIAL TOOLS & TEST EQUIPMENT:

Standard Field Engineer Tool Kit.
Anti-static Wrist Strap & Anti-static Work Mat.

FURNISHED MATERIALS:

FMI kit 2249157 consists of the following material and documentation:

Part Number	Description	Qty
2249157FMI	FMI 25258 Document	1
2151622	2-Way Video Splitter	1

Attention

-----SPECIAL RETURN-----

THE VIDEO SPLITTER THAT IS REPLACED IN THIS FMI NEEDS TO BE RETURNED TO THE VENDOR(BLACK BOX). PLEASE REUSE THE BOX THAT WAS SHIPPED WITH THE NEW VIDEO SPLITTER TO PACKAGE AND RETURN THE OLD VIDEO SPLITTER. A PRE-PRINTED FEDEX FORM HAS BEEN INCLUDED WITH THE KIT(by the vendor) TO BE ATTACHED TO THE BOX AND DROPPED OFF AT THE NEAREST FEDEX STATION.

PROCEDURE:

Introduction

Section 1.0 General Precautions

Section 2.0 Video Splitter Replacement

Section 3.0 System Check-out


Section 4.0 FMI Completion

SECTION 1.0 General Precautions

In the General System Service Manual(2152918-100), turn to Chapter 8, Paragraph 2.1.2.6 and read and become familiar with the General Precautions for handling and working on the Octane computer Hardware.

SECTION 2.0 2-Way Video Splitter Replacement

1. Bring down applications and turn console power off.
2. Locate the placement of the 2-way splitter. It should be near the console.
3. Remove console front panel. (If applicable)
4. Unplug the video splitter power source.
5. Locate video splitter on or near the console and unplug all video coax connectors attached to the video splitter. Note which video cables were input and which ones were output.

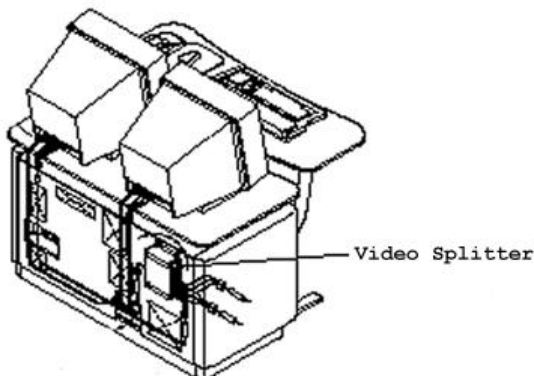

NOTE: This picture depicts only one possible location.

6. Remove the old video splitter and replace with new video splitter.
7. Reconnect all video cables paying close attention which ones were input and which ones were output.
8. Plug one end of the power cord into the back of the splitter and the other end into the console. (If applicable, because the video splitter could be plugged in elsewhere.)
9. Replace front cover of the console and power back up.

SECTION 3.0 System Checkout

Verify that all monitors connected to the video splitter are displaying images.

SECTION 4.0 FMI Completion

1. Repackage video splitter and attach the FED EX label included and drop off at your nearest FED EX pick up station. **(THIS IS A MUST FOR COMPLETION OF THE FMI!)**
2. Fill out and debrief the FMI Completion Dispatch Report in accordance with your local procedure.

g

GE Medical Systems

CT SERVICE ENGINEERING, NB-918

P.O. BOX 414

MILWAUKEE WISCONSIN 53201-0414 U.S.A.

Fax: (414) 785-5874