

**ADMINISTRATIVE INFORMATION****FMI No: 25298****FMI Title: LightSpeed Ultra Regulatory Gantry Tilt Covers Upgrade**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☒ Batch with	25300, 25307
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked: 228117****FMI Instructions Part No.: 2324555-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 2529825298

FMI Title: LightSpeed Ultra Regulatory Gantry Tilt Covers Upgrade

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments _____

Effectivity

CT LightSpeed Ultra Systems

Gantry Model Number: 2281177

Serial Numbers: 279519CN0 through 286390CN0.

Purpose

The purpose of this FMI is to install additional covers to the base of the gantry. When tilted and rotating the rotate frame of the gantry is exposed and presents a physical hazard. Safety and Regulatory Engineering has determined that the addition of these covers is necessary to protect the operator.

The gantry covers added by this FMI will need to be removed before performing a Detailed Calibration and then reinstalled at the end of the process. The Detailed Calibration software (see FMI25300) presents a warning message indicating the need for cover removal via a pop up window.

On some systems shipped during Fiscal Weeks 42, 2001 through 47, 2001 the gantry fan top cover limit switches were set too low (50 mm). Be sure to adjust the limit switch height, if needed, as indicated in **Gantry Fan Top Cover Limit Switch Adjustment Procedure** section.

A few systems may have shipped with incorrect gantry thermostat settings of 28, 1. Please be sure to verify your thermostat settings are 26, 2.

Prerequisite FMI's

- 25300

Related FMI's

- None

Tools Required

Qty	Description	Part #
	None	

People Power

Primary FE for .5 hour labor and 1 hour of travel.

Furnished Materials (2325443)

Part #	Description	Qty
2324555-100	FMI Document	1
2294827	Left Side Panel	1
2294827-2	Right Side Panel	1
46-328417P12	Screw 6mm X 12mm	8
46-328431P2	Washer, Plain 6.4mm X 18mm	8
46-328432P3	Washer, Lock 6.1mm X 11.8mm	8

Procedure

Right Side Panel Installation Procedure (2294827-2)

1. Remove gantry side covers.
2. Turn off axial drive enable.
3. Tilt gantry forward to around +18.5 degrees.
4. Referring to Figure 1, carefully slide new side panel (2294827-2) into position from the rear of the gantry.
5. Tilt gantry backward to around -5.0 degrees.
6. Disconnect right lower gantry front cover latch.

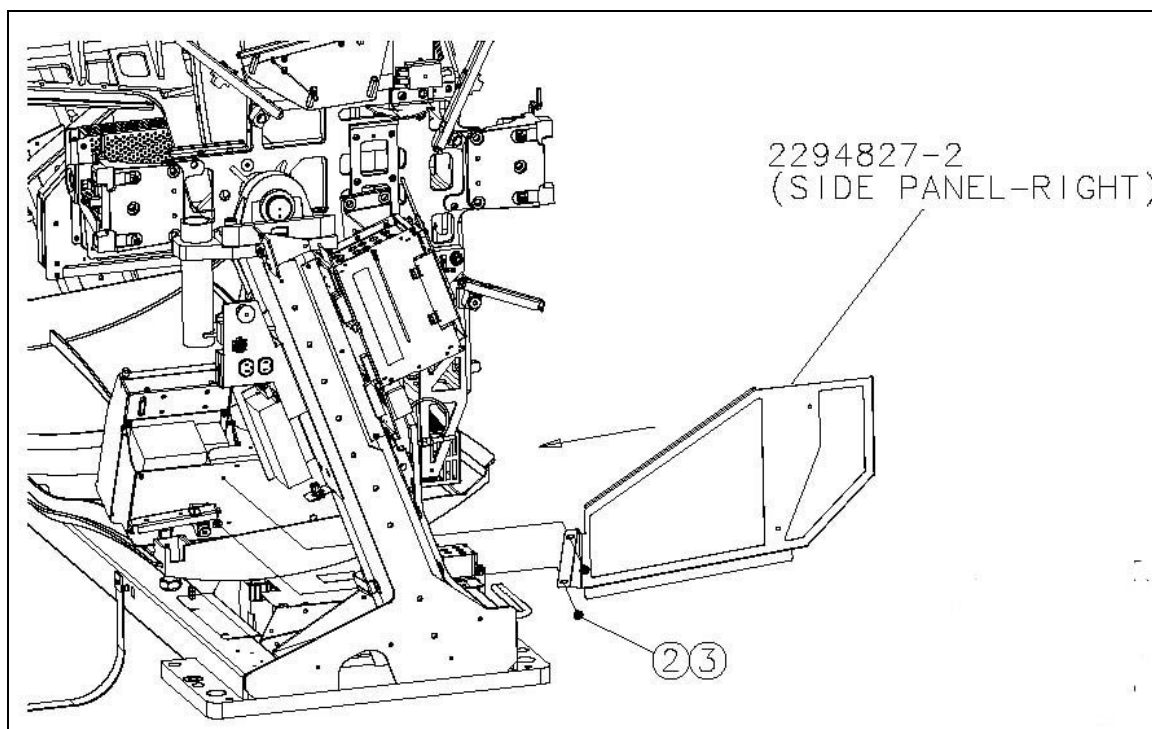


Figure 1. Right Side Panel Placement

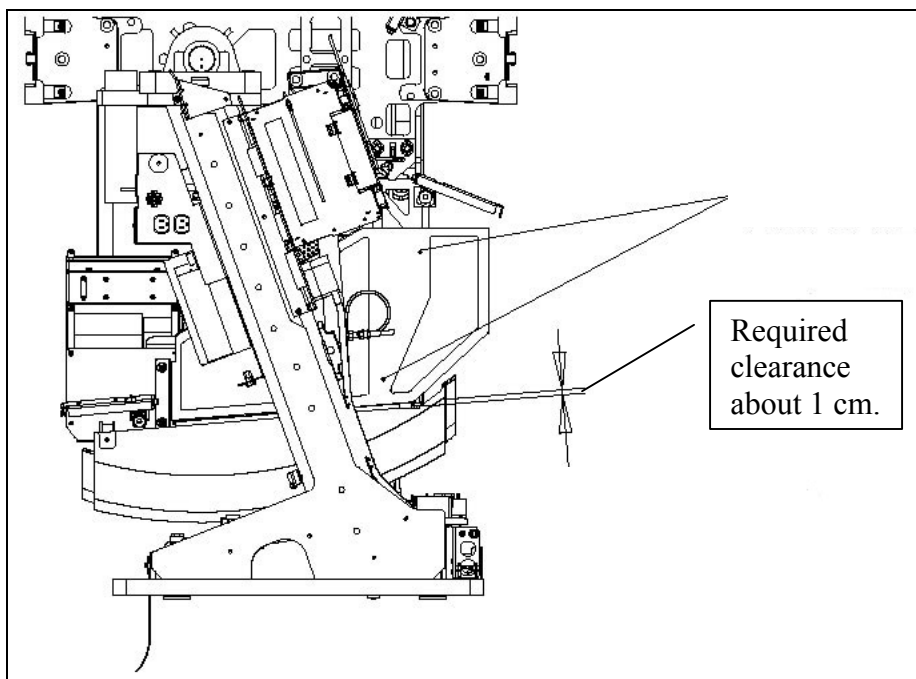


Figure 2. Right Side Panel Alignment

7. Loosely fasten two front screws (46-328417P12), two lock washers (46-328432P3), and two washers (46-328431P2) to the new holes at the front edge of the panel (see Figure 2).
8. Loosely fasten two rear screws (46-328417P12), two lock washers (46-328432P3), and washers (46-328431P2) to the new holes at the rear of the panel (see Figure 2).
9. Assure that the clearance indicated in Figure 2 (1 cm) is met before final tightening of the four screws.
10. Tighten screws.
11. Reconnect right lower gantry front cover latch.

Left Side Panel Installation Procedure (2294827)

1. Verify axial drive is turned off.
2. Tilt gantry to 0 degrees.

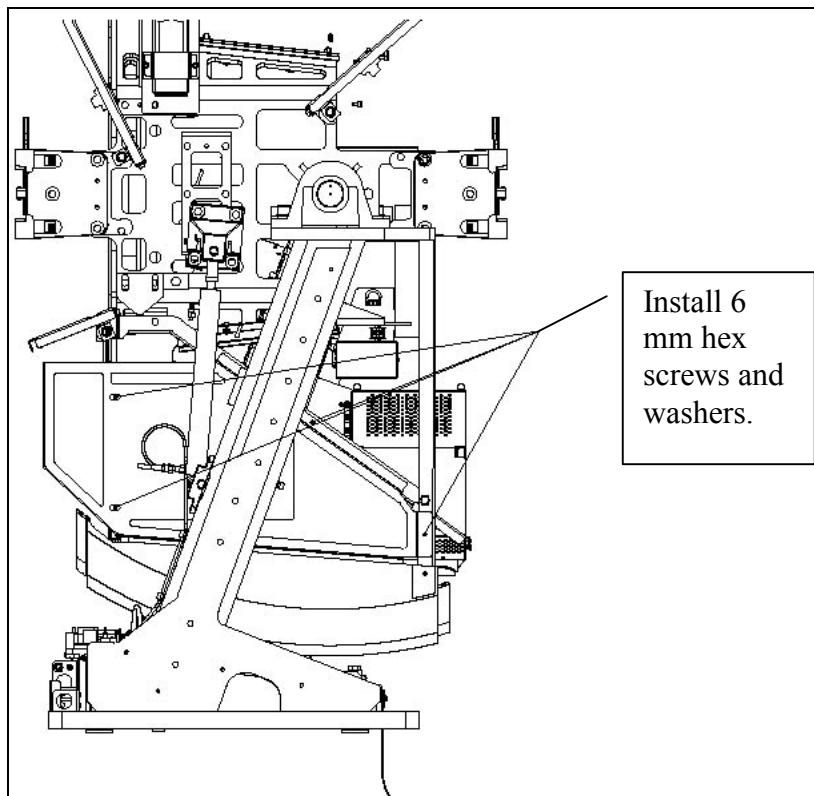


Figure 3. Left Side Panel Installation

3. As shown in Figure 3, install Left Side Panel (2294827) using four screws (46-328417P11), four lock washers (46-328432P3), and four washers (46-328431P2) supplied in FMI kit.

Gantry Fan Top Cover Limit Switch Adjustment Procedure

1. Verify and shut off gantry 120 VAC.
2. Remove top covers.
3. Per Figure 4, verify the limit switch height is set to 55 mm. If not, adjust the switch plunger height to this distance.

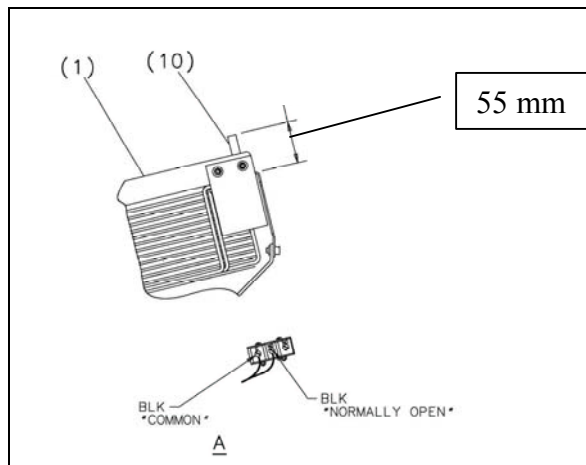


Figure 4. Limit Switch Height Adjustment

4. Place top covers onto gantry and verify limit switches engage. Referring to Figure 4, adjust as need to obtain positive actuation of the switches.

Finish Up

1. Turn all power back on.
2. Verify the gantry thermostat settings are "26" and "2."
3. Install all covers.
4. Inspect all gantry cover attachments for proper fit and latching.
5. Reset system.
6. Tilt gantry to 0 degrees and verify that the covers do not scrap or interfere with stationary parts of the gantry.
7. Tilt gantry to +30 and -30 degrees and verify that the covers do not scrape or interfere with stationary parts of the gantry.
8. Verify scanner is ready for customer use.
9. Recycle left over materials.