

**ADMINISTRATIVE INFORMATION****FMI No: 25310A****FMI Title: LS2002 M3 SDAS Software FMI**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety ☐ Mandatory ☒ Mandatory on Request	GEMS-AM: GEMS-E: GEMS-A:	Issue Date: Close Date:

☐ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole. ☐ Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole. ☐ Parts Ordering: GEMS - Americas, contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator ☐ No Parts are required to complete this FMI.
☐ Batch with	To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2353923-100****Do NOT leave unsecured on site.**



Completion Certification Report

Mail To

GEMS - Americas

GE Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe

GEMS - E / ESC
FMI Admin. 322
BP34
F 78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25310A FMI Title: LS2002 M3 SDAS Software FMI

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code:

- 1.) Modification Complete or previously modified.
- 2.) FMI not required per FMI instruction.
- 3.) Not modified. Unit in storage. FMI kit returned to stock.
- 4.) Modification refused or equipment altered. Customer signature and title required.
- 5.) Unit location unknown or unit not in area. Product Locator notified.

**Customer Refusal or Equipment Altered
- Modification Pending.**

Installer
Comments

Effectivity

1 "First Wave"

For the "first wave" of implementation, LightSpeed Plus (SDAS) systems with the following Gantry system IDs are affected:

System ID	Gantry Model #	Gantry Serial Number
910212091	—	—
920219019	—	—
HC2707CT01	2260619	00000269220CN0
083021592012120	2260619	00000270628CN1
—	2180482	0000050572CN5

2 "Second Wave"

For the "second wave" of implementation, LightSpeed Plus (SDAS) systems with the following Gantry system IDs are affected:

System ID	Gantry Model #	Gantry Serial Number
304327LSP	2260619	269038CN6
617754CT01	2260619	283114CN7
M5484926	2260619	00000267917CN3
M4189241	2260619	00000268056CN9
M4047048	2260619	00000272235CN3
M4548524	2260619	00000283009CN9
A5183747	2260619	00000265771CN6
A5333013	2260619	00000265016CN6
QY0002	2260619	00000262387CN4
QY0008	2260619	00000262929CN3
QY0029	2260619	00000278656CN4
QY0044	2260619	00000279055CN8
850210465	2260619	00000270418CN7
850210480	2260619	00000282187CN4
954721LS1	2260619-3	843216YM3
10331CTS49	2260619	00000262282CN7
00131CTS03	2260619	00000265295CN6
70282CTS02	2260619	00000265880CN5
00441CTS01	2260619	00000265942CN3
70638CTS02	2260619	00000265943CN1
05451CTS01	2260619	00000269281CN2
00298CTS05	2260619	00000273431CN7
10079CTM51	2260619	00000276416CN5
10140CTM54	2260619	00000276737CN4

System ID	Gantry Model #	Gantry Serial Number
70000CTM05	2260619	00000277958CN5
00362CTS01	2260619	00000280122CN3
573348LCT	2260619	276224CN3
847360QXI	2260619	282185CN8
8987021796	2260619	00000264895CN4
910213050	2260619	00000269191CN3
910213046	2260619	00000261444CN4
910213049	2260619	00000266834CN1
082421040053	2260619	00000264393CN0
850210489	2260619	00000280124CN9
0834210012	2260619	00000261759CN5
0856210051	2260619	00000268647CN5
0847210009	2260619	00000263807CN0
853210112	2260619	00000273900CN1
853210113	2260619	00000279528CN4
853210115	2260619	00000279517CN7
580220LSP	2260619	00000278447CN8
920219019	2260619	00000262167CN0
910212090	2260619	00000277980CN9
10097CTM55	—	—
910212077	—	—
910525010	—	—
910213053	—	—
910213052	—	—
082421030103	—	—

Purpose

The purpose of this FMI is to support customers who are level 4 escalations and customers who were promised Cardia Burst and Smart Step.

Furnished Materials

FMI 25310 Kit (Part Number 2354013) includes the following items:

Qty	Description	Part #
1	FMI Document	2353923-100
1	LS2002 SW Set (SDAS) (see list, below)	2350063

The LS2002 SW Set (SDAS) includes the following items:

Qty	Description	Part #
1	LS2002 H2 M3 OS and Apps (SDAS)	2348410
1	LS2002 Software Installation Manual	2348269-100
1	Application Tips and Workarounds	2243321-100
1	2.0 LightSpeed Plus Tech Ref Manual	2281207-100
1	2.0 LightSpeed Plus Operator's Manual	2281208-100
1	Serial #, Model #	46-302200P7
1	440 mA Scan Option	2345411

The LS2002 H2 M3 Software Kit (P/N 2348410) includes the following items:

Qty	Description	Part #
1	Operating System Disk 1	2305087-2
1	Operating System Disk 2	2305088
1	Operating System Disk 3	2305089
1	Applications Software	2348409
1	LS2002 M3 Patch CD	2349518


Required Tools

Qty	Description	Part #
1	LightSpeed Plus (SDAS) Class C	2344426
1	LightSpeed Plus (SDAS) Class M	2344427

People Power

4 hours labor and 1 hour of travel.

Steering Guide

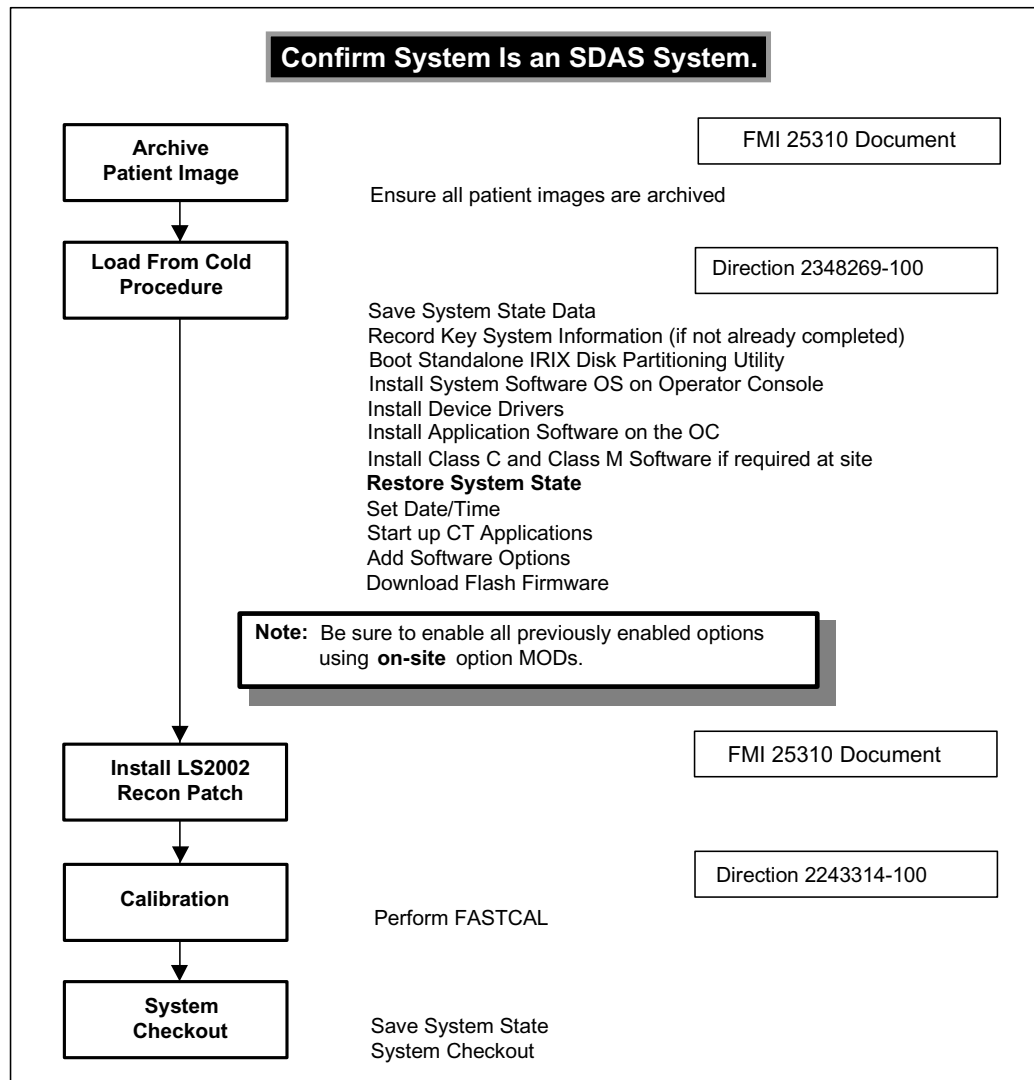

Procedure

1 Archive Patient Images

Verify that the customer has archived all patient data.

2 Software Load Procedures (LFC)

Perform the Load From Cold (LFC) procedure as detailed in Chapter 1 of *LS2002 Software Installation Procedures* (2348269-100).



NOTICE

It is required that a “Restore STATE” be performed following the load of Class C and Class M software. To facilitate this, it is recommended that the Class C and Class M be loaded following the Applications Load (Section 3.6).

3 Installing the LS2002 M3 Recon Patch

3.1 Installation

- 1.) Shut down applications by:
Service -> Utilities -> Application Shutdown
- 2.) Open a UNIX shell and become root by:
type: `su -`
- 3.) Enter the root password.
- 4.) Insert the LS2002 M3 Recon Patch CD in the CD-ROM Drive:
 - a.) type: `mountCD /CDROM`
 - b.) type: `/CDROM/install_ls2002_m3_patch`The install messages will appear:
.....
.....
- 5.) When the installation is complete:
 - a.) type: `unmountCD`
 - b.) Reboot the system

3.2 Verification of the Software Patch Installation

To verify Patch is installed in OC:

type: `showprods | grep ls2002_patch_b`

The output will look similar to the following:

```
ls2002_patch_b 07/10/2002 Patch for ls2002 M3
ls2002_patch_b.img 07/10/2002 Patch files for recon
ls2002_patch_b .img.ocroot 07/10/2002 Files that go onto ocroot
```

4 Calibration

Perform FASTCAL by selecting Daily Prep on the bottom of the left monitor, then selecting FASTCAL from the displayed selections.



NOTICE You may get a pop-up indicating: “System detects an alignment or change was performed on the System”. If this happens, press the Ok button and continue with the FASTCAL.

5 System Checkout and Save State

- 1.) Refer to the Software Installation Procedures Manual, and perform a Save System State. This will save the new calibrations as well as the customer and service protocols.
- 2.) Refer to the Software Installation Procedures Manual, and perform the system performance checks as described in the manual. Ensure all networking and filming processes are still functional.
- 3.) Complete and return Red Installation Product Locator Cards.

6 Appendix: Uninstalling the LS2002 M3 Recon Patch

- 1.) Shut down applications by:
Service -> Utilities -> Application Shutdown
- 2.) Open a UNIX shell and become root by:
type: `su -`
Enter the root password.
- 3.) Insert the LS2002 M3 Recon Patch CD in the CDROM Drive:
 - a.) type: `mountCD /CDROM`
 - b.) type: `/CDROM/uninstall_ls2002_m3_patch`The uninstall messages will appear:
.....
.....
- 4.) When the installation is complete:
 - a.) type: `cd`
 - b.) type: `unmountCD`
 - c.) Reboot the system