

**ADMINISTRATIVE INFORMATION****FMI No: 25360****FMI Title: Helios Collimator Xray Scatter Filler Bar**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support		
GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked: See "Effectivity" Section****FMI Instructions Part No.: 5133567-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25360

FMI Title: Helios Collimator Xray Scatter Filler Bar

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments

Effectivity

All CT and PET/CT system build prior to 1/19/2005 with following gantry part numbers:

For Hino 221478-3 collimator impacted gantry p/n's:

- 2341799
- 2341799-2
- 2341799-3
- 2260619-3
- 2281177-3
- 2281177-4
- 2281177-10
- 2281177-12

For Hino 221478-4 collimator impacted gantry p/n's:

- 2378661
- 2378661-3
- 2378661-10
- 2401904
- 5120327
- 2357739-7
- 2383550-2
- 2383550-10
- 2383549-2
- 2383549-10
- 5119187
- 5127644-10
- 2401904

For GEMSAm 221478-3 collimator impacted gantry p/n's:

- 2180482
- 2294396
- 2281177
- 2281177-2
- 2281177-7
- 2339985
- 2339985-2
- 2339985-4
- 2380553

For GEMSAm 221478-4 collimator impacted gantry p/n's:

- 2357739
- 2357739-3
- 2357739-4
- 2357739-6
- 2374681
- 2374681-2

Purpose

The purpose of this FMI is to install a small filler bar (p/n 5134695) on the collimator frame to reduce the risk of potential X-ray containment issues at the seam between the filter assembly and the collimator frame.


Prerequisite FMIs

None

Related FMIs

None

Furnished Materials

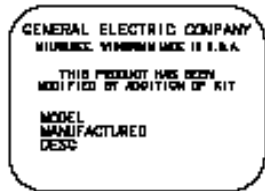Qty	Description	Part #	Comments
1	FMI Document	5133567-100	
1	M4 x 0.7 bottoming-style hand tap	5135684	
1	Rating Plate - Sticker	46-302200P12	
1	Notched Filler Bar Assembly	5134695	
1	FMI 25360 Installation-Complete Letter	5141853-147	FE to deliver to Radiology Manager
1	Important Radiation Safety Notice	5141854-147	FE to deliver to Radiology Manager

Tool Required

No additional tools; standard FE tool kit required

Personnel Requirements

1/2 hour labor

1 hour travel

Procedure

1. Remove side, top, and front gantry covers.
2. Turn off the gantry axial enable, HVDC enable and 120 Vac switches.
3. Lockout all power.
4. Rotate the gantry so that the collimator is at the 6 o'clock position.
5. As a precautionary measure, cover the collimator output port with a towel to prevent debris from getting onto the port.
6. Remove the (2) existing 4mm screws from the front of the collimator filter mounting plate. Discard the screws.



Figure 1 - Example of screw being removed



Figure 2 - Example of screw being removed

7. Using the provided hand tap (p/n 5135684), clean the mounting holes of LocTite by gently turning the tap in each mounting hole until it lightly bottoms out.
To avoid cross-threading or breaking the tap, do not force.
8. With the notch facing downwards, attach the filler bar assembly (5134695) to the collimator by starting to install the two filler bar assembly screws into the empty screw holes on the collimator.

NOTE: The filler bar only fits in one direction – notch downward.

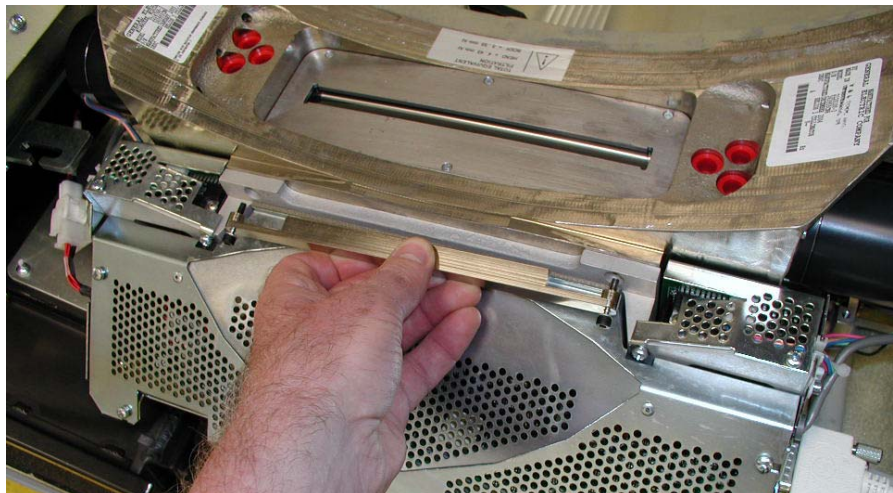


Figure 3 - Installing filler bar

9. Torque the two M4 screws on the notched filler bar assembly to 2.3 Nm (20.4 lb-in) using the appropriate torque wrench.

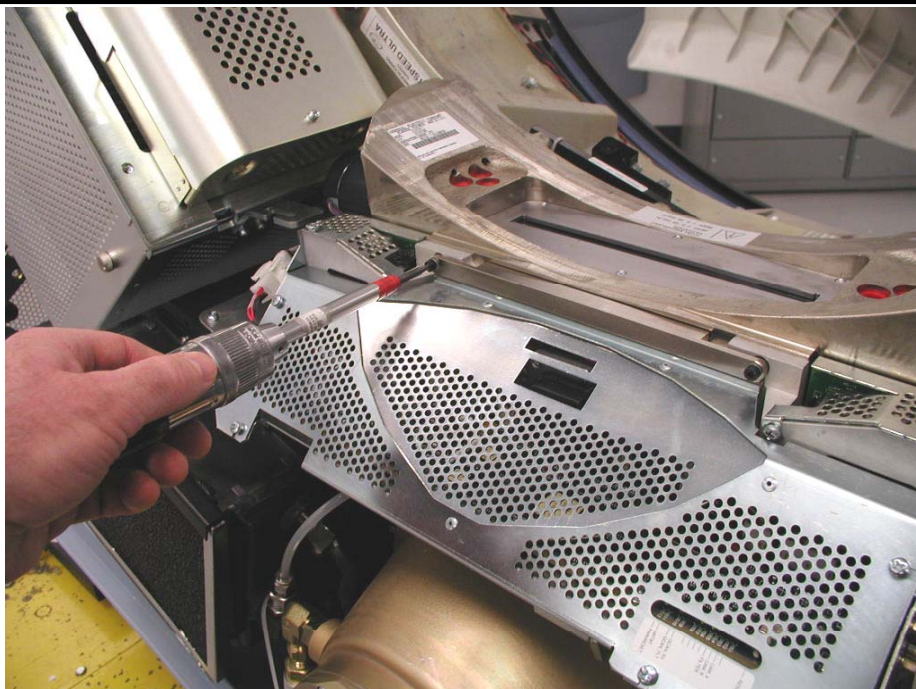


Figure 4 - Torque the M4 Screw



Figure 5 - Torque the M4 Screw

10. Place the rating plate sticker from this FMI kit on the collimator frame as shown in the picture.

Do not cover up any text on the existing labels. Folding the sticker over an edge is acceptable, although this is not shown in Figure 6.



Figure 6 - Rating Plate Sticker Placement

11. Installation is complete.

Finishing Up

Upon completion of this FMI, deliver the letters included in this FMI kit to the radiology manager or department supervisor:

- 5141853-147, FMI 25360 Installation / corrective action-complete letter
- 5141854-147, Important Radiation Safety Notice

Commonly used "Hazard Icons" icons:

