

**ADMINISTRATIVE INFORMATION****FMI No: 25368****FMI Title: VCT PDU Update**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS - AM:	Issue Date: Feb. 3rd, 2006
☒ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date: Sep. 1st, 2006

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 5144784-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas

GEMS - Europe

G.E. Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe
FMI Admin. 322
BP34
F78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: FMI 25368

FMI Title: VCT PDU Update

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments

Effectivity

Existing VCT systems:

Froedtert	414805CT4
CCN	M2322223
JFK	732321VCT
Keio University Hospital	XG0009
Hong Kong Sanatorium	0828210048
Univ. of Washington	206598VCT2
Charite	CC0618CT07
AsiaMedic	0847210010
Toki Municipal Hospital	QS0001
Handai Fuzoku (Osaka)	QX0008
Dortmund	HC0074CT04
Huntsville Medical Center, AL	256MAINVCT
Cardiology Associates, AL	251340VCT
West Suburban, IL	630972VCT
CHUV Lausanne	260010CT15
RIKSHOSPITALET HF	N01016CT06
University of Michigan	734936VCT
UNIVERSITA' DEGLI STUDI DI ROMA TOR VERGATA, Italy	A5333029
AMI	609652VCT
Miami Baptist, FL	305596VCT
Boston Med Ctr, MA	617414VCT1
North Shore Univ., NY	516562VCT
Bronson Methodist Hospital, MI	269341VCT
Mobile Infirmary	251435CT3
IMorristown, NJ / Mountainside	973429VCT
East River Medical imaging, NY	212288CT64
St. Mary's Hospital, MI	989497CT
Methodist Hospital, TX	713441MCT9
HOPITAL BEAUJON	M4016557
LHSC - UNIVERSITY HOSPITAL	519685VCT

Showa Univ Fujigaoka	QV0026
IRYOU HOUJIN SHADAN HOKUTO HOSPITAL	QV0030
HOPITAL NECKER ENFANTS MALADES	M4033088
KEN KOUSEI NOUKYOU RENGOU HIROSHIMA SOUGOU HOSPITA	QS0002
HOPITAL EUROPEEN GEORGES POMPIDOU	M5460685
KENKOU HOKEN ISAHAYA SOUGOU HOSPITAL	QS0005
UNIV OF WASHINGTON MED CTR	206598VCT1
MC Morice (Massy)	M9227209
Metropolitan Health Service (Sir Charles Gardiner)	0910216045
Huddinge (Karolinska)	A007CT01
KOBE SHIRITSU CHUO SHIMIN HOSPITAL	QS0007
INOVA MT VERNON HOSPITAL	703664LS64
PARK RIDGE HOSPITAL	828681CT
EVANSTON NORTHWESTERN HEALTHCARE CORP	847570VCT
Cardiff	00243CTS08
ROSWELL PARK CANCER INSTITUTE	716845VCT
BAYLOR HEART AND VASCULAR	214820VCT
ED CENTER--Service Training + a console	414541VCT2
Lifecare Diagnostic, Malaysia	0834210025
BETH ISRAEL DEACONESS MEDICAL CENTER	617754VCT7
Zurich	900010CT04
ST JOSEPHS HEALTH CENTRE	519646VCT
AKRON GENERAL HEALTH SYSTEM	330996VCT
Military Warsaw	600042CT02
KAISER HOSPITAL	408972VCT1
QUEST IMAGING MED ASSOC INC	661633VCT
Gottingen (Grabe)	UC1220CT10
NORTH MISSISSIPPI MEDICAL CENT	662CROSSCT
IRYOU HOUJIN KOUJINKAI HOSHIGAURA HOSPITAL	QS0006
Monaco	X751958013
Polyclinic di Monza	A5401345
BOSTON CITY HOSPITAL	617414VCT2
PHSA(QC)/QUEEN ELIZABETH HOSP	902894VCT
ST BERNARDS MEDICAL CTR	870972VCT
CANTERBURY DISTRICT HEALTH BOARD	0920 21 9032

ESKILSTUNA	D007CT03
"DRS LAKE, SMIT & PARTNERS"	ZA1313CT02
Shinkoga Hospital	QS0008
ST LUKES MEDICAL CENTER	414649VCT5
OAKWOOD HOSPITAL, Dearborn, MI	313593VCT
ST JOSEPHS HEALTH CENTRE	416534VCT1
Chiba Rosai	QS0003
M.A.R. LTD.	IL1086CT03
SCM RADIOLOGIE IMAGERIE MEDICALE	M1322107
HCA INTERNATIONAL LTD	00888CTS01
HARALDSPLASS DIAKONALE SYKEHUS	NO1146CT01
AZ MARIA MIDDELARES	110014CT04
STRAHLEN INSTITUT CDT GMBH	PC6066CT05
CHRISTUS ST MICHAEL HLTH SYSTE	903614VCT1
Matsuzaka Chuo Shimin Hp	QS0009
BUFFALO GENERAL HOSPITAL	716859VCT
CMC DE BLIGNY	M9024408
Penang Adventist Hospital	0834210026

Purpose

The purpose of this FMI is to update early VCT PDUs to "Version B"; approximately 80 early systems are affected. The PDU was redesigned during the VCT program to provide added protection in the event that an input phase would drop. The new design will turn off all output power to the system if one of the main input phases drops out.

Prerequisite FMIs

None

Related FMIs

None

Furnished Materials

FMI kit (5144827)

Qty	Description	Part #
1	FMI Document	514478-100
1	NGPDU-4 Assembly	2326492-4

Tool Required

Qty	Description
1	Standard field service kit


Personnel Requirements

"2" hours labor and "1" hours of travel

Procedure

Preparation, PDU Inspection

1. Verify that the existing PDU is the correct version by looking for a sticker on top right portion of frame that refers to "version B". If the label is present then verify 2 green voltage monitor relays are present on the DIN rail, in the middle of the circuit breaker assembly.

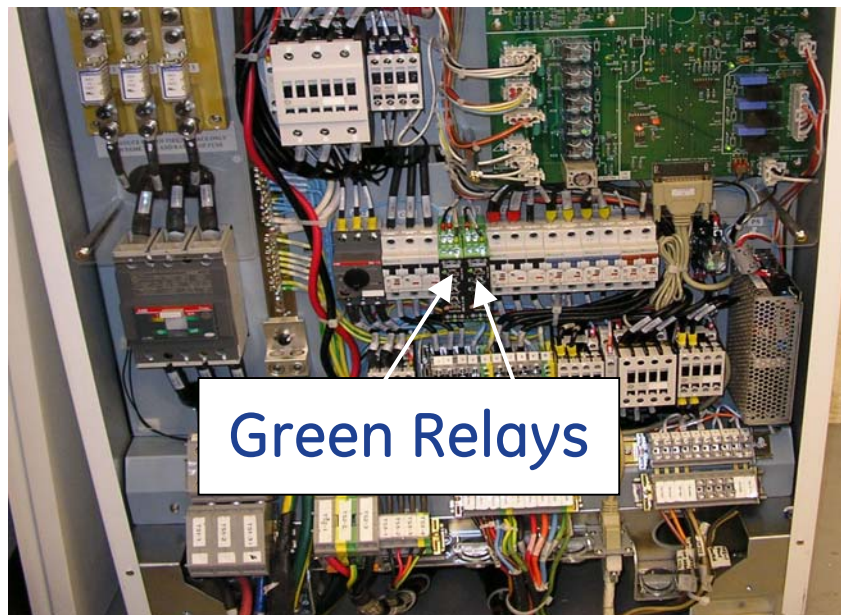

2. If the "Version B" label is present and the green relays are not, please contact Eric Stepanovich @ 262-312-7736.
3. If the existing PDU **does not have the label with "version B", the PDU will need to be replaced with the kit provided.** Proceed to the PDU replacement section below.
4. If the existing PDU **has the label with "Version B" on it and the 2 green voltage monitor relays are present this is the correct PDU** stop here. The FMI does not need to be performed. Return unused FMI kit to: Call 262-754-0424 with shipping details.

HM Product Solutions, LTD.

3225 Gateway Road, Suite 250

Brookfield, WI 53045

262-754-0424

PDU Replacement

1. Summary:
 - Verify PDU version.
 - Remove Old PDU – Install New PDU
 - PDU Initialization and configuration
 - Finalization
2. For those sites that require PDU replacement, please refer to connection configurations detailed in the PDU Cable Connections & Configuration, in Section 7.0 of LightSpeed 7.X Installation Manual.

The documentation can be found in the LS 7.X Service Methods CD (5116423)

- Click on Install (Dir 5116411-100)
 - Click on LightSpeed 7.X Installation Manual
 - Select Book 1 – Mechanical Installation
 - Go to Chapter 2 - Power, Ground & Interconnect Cables
 - Then to Section 7.0 – PDU Cable Connections & Configuration
 - Follow accordingly
3. To complete the PDU installation refer to the LS 7.X Service Methods CD (5116423)
 - Click on Install (Dir 5116411-100)
 - Click on LightSpeed 7.X Installation Manual
 - Select Book 2 – Electrical Calibration, Integration and Testing
 - Go to Chapter 5 - Electrical Integration and Safety Verifications
 - Then to Section 1.4 – Initial PDU Configuration
 - Follow accordingly
 4. Go to Finalization

Finalization

1. Perform a test scan using the 20cm QA phantom to verify system operation.
2. **After the PDU has been replaced, it must be returned to the FMI supplier at the following address, See attached shipment documentation:** Call 262-754-0424 with shipping details.

Note: It is very important that the PDU is returned as quickly as possible and in working order. Please do not take parts off of the PDU that is being returned, it will be reworked and shipped to the next site.

HM Product Solutions, LTD.
3225 Gateway Road, Suite 250
Brookfield, WI 53045
262-754-0424

3. Close FMI