

**ADMINISTRATIVE INFORMATION****FMI No: 25372****FMI Title: VCT M3A (05MW31.6) Software update**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date: 9/23/2005
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date: 11/25/2005

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 5144790-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas

GEMS - Europe

G.E. Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe
FMI Admin. 322
BP34
F78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25372

FMI Title: VCT M3A Software update

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments

Effectivity

The following CT Equipment is affected:

MODEL	SERIAL	GON	Site		SYSTEM	POLE
5124069	00000037107CN5		EVANSTON	IL	847570VCT	AMER
5124069	00000359252CN4		MILWAUKEE	WI	414805CT4	AMER
5124069	00000359394CN4	532604	METUCHEN	NJ	732321VCT	AMER
5124069	00000362689CN2	544256	BUFFALO	NY	716845VCT	AMER
5124069	00000365378CN9	553702	BUFFALO	NY	716859VCT	AMER
5124069	00000366067CN7	537379	SAN DIEGO	CA	858966CT1	AMER
5124069	00000366278CN0	562470	SEATTLE	WA	206598VCT2	AMER
5124069	00000367275CN5	544545				AMER
5124069	00000367835CN6	559698				AMER
5124069	00000367853CN9	153882	ANN ARBOR	MI	734936VCT	AMER
5124069	00000368182CN2	151000	NORTHFIELD	NJ	609652VCT64	AMER
5124069	00000368454CN5	545670	MIAMI	FL	305596VCT	AMER
5124069	00000368496CN6	554278	MANHASSET	NY	516562VCT	AMER
5124069	00000368498CN2	568013				AMER
5124069	00000368597CN1	175442	BOSTON	MA	617414VCT1	AMER
5124069	00000368785CN2	551660	KALAMAZOO	MI	269341VCT	AMER
5124069	00000368787CN8	547010	MONTCLAIR	NJ	973971VCT	AMER
5124069	00000368796CN9	571473	NEW YORK	NY	212288CT64	AMER
5124069	00000368864CN5	152707				AMER
5124069	00000368995CN7	540176	HOUSTON	TX	713441MCT9	AMER
5124069	00000369006CN2	192133				AMER
5124069	00000369119CN3	652594				AMER
5124069	00000369976CN6	570902				AMER
5124069	00000370180CN2	579137	BIRMINGHAM	AL	205783VCT	AMER
5124069	00000370181CN0	522420				AMER
5124069	00000370189CN3	320212	ONTARIO	CA	519685VCT	AMER
5124069	00000370257CN8	562467				AMER
5124069	00000370521CN7	376596	FLETCHER	NC	828681CT	AMER
5124069	00000370890CN6	553146	ALEXANDRIA	VA	703664LS64	AMER
5124069	00000371017CN5	553857	EVANSTON	IL	847570VCT	AMER
5124069	00000371094CN4	50513				AMER

MODEL	SERIAL	GON	Site		SYSTEM	POLE
5124069	00000371101CN7	566462	DALLAS	TX	214820VCT	AMER
5124069	00000371240CN3	548901	AKRON	OH	330996VCT	AMER
5124069	00000371401CN1	438244				AMER
5124069	00000371402CN9	570620	LONDON ON CAN		519646VCT	AMER
5124069	00000371436CN7	528381	BOSTON	MA	617754VCT7	AMER
5124069	00000371592CN7	156300	CHARLOTTETOWN PE	CA	902894VCT	AMER
5124069	00000371636CN2	573354	SAN JOSE	CA	408972VCT1	AMER
5124069	00000371671CN9	566607	TUPELO	MS	662CROSSCT	AMER
5124069	00000371809CN5	542720	BOSTON	MA	617414VCT2	AMER
5124069	00000371957CN2	484353	JONESBORO	AR	870972VCT	AMER
5124069	00000372432CN5	50606				AMER
5124069	00000372502CN5	571020	DEARBORN	MI	313593VCT	AMER
5124069	00000372516CN5	531554	TORONTO ON CAN		416534VCT1	AMER
5124069	00000372618CN9	573104	WAUWATOSA	WI	414649VCT5	AMER
5124069	00000372633CN8	556231				AMER
5124069	00000372677CN5	559339	TEXARKANA	TX	903614VCT1	AMER
5124069	00000372861CN5	574412	KNOXVILLE	TN	865541CT2	AMER
5124069	00000372938CN1	570905	ATLANTA	GA	404712VCT	AMER
5124069	00000372944CN9	568415	GUELPH ON CAN		519822VCT	AMER
5124069	00000372974CN6	643704	CINCINNATI	OH	513585CT2	AMER
5124069	00000372976CN1	574577	DANVILLE	PA	814231VCT	AMER
5124069	00000373040CN5	572830	HATTIESBURG	MS	601288VCT1	AMER
5124069	00000373053CN8	643782	KALAMAZOO	MI	269381VCT	AMER
5124069	00000373054CN6	560713	SUGARLAND	TX	281274VCT	AMER
5124069	00000373199CN9	558929				AMER
5124069	00000373266CN6	560749	PHILADEPHIA	PA	215456VCT	AMER
5124069	00000373272CN4	533369	WALKERTON ON CAN		519881VCT	AMER
5124069	00000373388CN8	628600	KNOXVILLE	TN	865545VCT	AMER
5124069	00000373396CN1	578408	WINSTON SALEM	NC	336713EDVCT	AMER
5124069	00000373648CN5	579216	LAFAYETTE	LA	337261UVCT	AMER
5124069	00000373649CN3	691527				AMER
5124069	00000373652CN7	774207	MONTREAL PQ CAN		514252VCT	AMER
5124069	00000373654CN3	578026	ASHLAND	KY	606833VCT	AMER
5124069	00000373659CN2	561682	PORTLAND	OR	503413AW	AMER
5124069	00000373701CN2	571809				AMER

MODEL	SERIAL	GON	Site		SYSTEM	POLE
5124069	00000373702CN0	562529	AURORA	IL	630801VCT1	AMER
5124069	00000373705CN3	571505	RICHMOND	VA	703689VCT	AMER
5124069	00000373719CN4	568859	HARTFORD	CT	860545VCT	AMER
5124069	00000373861CN4	568437	LOS ANGELES	CA	310278VCT	AMER
5124069	00000373911CN7	578182	MISSOULA	MT	406543VCT	AMER
5124069	00000373912CN5	575847				AMER
5124069	00000374169CN1	564427	LEBANON	NH	603650VCT	AMER
5124069	00000374532CN0	569939	MACON	GA	478633VCT	AMER
5124069	00000374544CN5	800257				AMER
5124069	00000374581CN7	567112	SAN FRANCISCO	CA	415353VCT1	AMER
5124069	00000374583CN3	570865				AMER
5124069	00000374589CN0	566661				AMER
5124069	00000374593CN2	580010	DAYTON	OH	937208VCT	AMER
5124069	00000374595CN7	532713				AMER
5124069	00000374611CN2	46163				AMER
5124069	00000374612CN0	582045	PHOENIX	AZ	602604VCT	AMER
5124069	00000374618CN7	706245				AMER
5124069	00000374656CN7	578072	DEARBORN	MI	217876VCT1	AMER
5124069	00000374657CN5	561742				AMER
5124069	00000374659CN1	574636	PEORIA	IL	309672VCT	AMER
5124069	00000374663CN3	578186				AMER
5124069	00000374664CN1	580766	QUAKERTOWN	PA	215538VCT	AMER
5124069	00000374670CN8	573316				AMER
5124069	00000374829CN0	583490				AMER
5124069	00000374831CN6	793547	LAKE ZURICH	IL	847842VCT	AMER
5124069	00000374837CN3	547020				AMER
5124069	00000375457CN9	561059	NEW YORK	NY	COR359874	AMER
5124069	00000375780CN4	568198	STRATFORD	ONT	519272VCT	AMER
5124069	00000375800CN0	572294	SAN ANTONIO	TX	210208MOBVCT	AMER
5124069	00000375871CN1	848302				AMER
5124069	00000375883CN6	572975				AMER
5124069	00000375884CN4	572975				AMER
5124069	00000375889CN3	887428				AMER
5124069	00000375938CN8	568407				AMER
5124069	00000375940CN4	165097				AMER

MODEL	SERIAL	GON	Site		SYSTEM	POLE
5124069	00000375945CN3	879260				AMER
5124069	00000376002CN2	566236				AMER
5124069	00000376009CN7	554316				AMER
5124069	00000376026CN1	938020				AMER
5124069	00000376058CN4	519662				AMER
5124069	00000376061CN8	569112				AMER
5124069	00000376080CN8	926434				AMER
5124069	00000376081CN6	575085				AMER
5124069	00000376086CN5	574091				AMER
5124069	00000376277CN0	574599				AMER
5124069	00000376289CN5	581984				AMER
5124069	00000376299CN4	571457				AMER
5124069	00000376301CN8	581079				AMER
5124069	00000376308CN3	943248				AMER
5124069	00000376343CN0	945382				AMER
5124069	00000376345CN5	944993				AMER
5124069	00000361977CN2	999999	TOKYO	TOKYO	XG0009	ASIA
5124069	00000364830CN0	529417				ASIA
5124069	00000366477CN8	543828	SINGAPORE 238868	SG	847210010	ASIA
5124069	00000366852CN2	131154	GIFU	GIFU	QS0001	ASIA
5124069	00000368381CN0	568385	BEIJING CITY		82421040151	ASIA
5124069	00000368402CN4	165189	BEIJING_CITY		82421040148	ASIA
5124069	00000369418CN9	577362				ASIA
5124069	00000370188CN5	560350	HOKKAIDO	HOKKAIDO	QV0030	ASIA
5124069	00000370201CN6	559047	TAIWAN		853210161	ASIA
5124069	00000370523CN3	547661	HIROSHIMA	HIROSHIMA	QS0002	ASIA
5124069	00000370578CN7	348248	NEDLANDS, WA	AU	910216045	ASIA
5124069	00000370627CN2	559544	NAGASAKI	NAGASAKI	QS0005	ASIA
5124069	00000370730CN4	50501	FUKUOKA	FUKUOKA	QS0008	ASIA
5124069	00000370751CN0	559554	HYOGO	HYOGO	QS0007	ASIA
5124069	00000370847CN6	569979				ASIA
5124069	00000371151CN2	562001	MALAYSIA 462000		834210025	ASIA
5124069	00000371287CN4	575604				ASIA
5124069	00000371729CN5	559504	HOKKAIDO	HOKKAIDO	QS0006	ASIA
5124069	00000371892CN1	555022	NZ		920219032	ASIA

MODEL	SERIAL	GON	Site		SYSTEM	POLE
5124069	00000372333CN5	563352	BEIJING CITY		82421040154	ASIA
5124069	00000372614CN8	576914				ASIA
5124069	00000372642CN9	563980	CHIBA	CHIBA	QS0003	ASIA
5124069	00000372663CN5	564725				ASIA
5124069	00000372833CN4	576677				ASIA
5124069	00000372838CN3	578881	MIE	MIE	QS0009	ASIA
5124069	00000372853CN2	553237	MALAYSIA		834210026	ASIA
5124069	00000372934CN0	571481	HEFEI CITY,ANHUI PROV	INCE	82421110055	ASIA
5124069	00000373050CN4	570649				ASIA
5124069	00000373270CN8	571088				ASIA
5124069	00000373335CN9	571548				ASIA
5124069	00000373862CN2	569777				ASIA
5124069	00000373920CN8	574493				ASIA
5124069	00000374540CN3	571486				ASIA
5124069	00000374628CN6	577785				ASIA
5124069	00000375772CN1	572863				ASIA
5124069	00000375783CN8	164131				ASIA
5124069	00000375804CN2	578614				ASIA
5124069	00000375838CN0	583243				ASIA
5124069	00000376027CN9	561465				ASIA
5124069	00000376089CN9	584372				ASIA
5124069	00000376284CN6	581506				ASIA
5124069	00000359588CN1	551108	SAINT DENIS		M2322223	EURO
5124069	00000366364CN8	990049	BERLIN 65		CC0618CT07	EURO
5124069	00000367119CN5	547193	DORTMUND 1		HC0074CT04	EURO
5124069	00000367354CN8	990041	ZUERICH	CH	900010CT04	EURO
5124069	00000367925CN5	537873	LAUSANNE		260010CT15	EURO
5124069	00000368022CN0	148890	OSLO	NO	NO1016CT06	EURO
5124069	00000368356CN2	565431	ROMA	IT	A5333029	EURO
5124069	00000370172CN9	548151	PARIS	FRANCE	M4033088	EURO
5124069	00000370363CN4	548510	PARIS		M5460681	EURO
5124069	00000370625CN6	572274	ANTONY		M9227209	EURO
5124069	00000370658CN7	560911				EURO
5124069	00000371087CN8	573696	CARDIFF	WA	00243CTS08	EURO
5124069	00000371288CN2	548058	GOETTINGEN		UC1220CT10	EURO

MODEL	SERIAL	GON	Site		SYSTEM	POLE
5124069	00000371432CN6	26564	WARSZAWA	PL	600042CT02	EURO
5124069	00000371814CN5	547369	MONACO		X751958013	EURO
5124069	00000372264CN2	560250	MILANO	IT	A5401345	EURO
5124069	00000372290CN7	569255	ESKILSTUNA		D007CT03	EURO
5124069	00000372311CN1	577265				EURO
5124069	00000372638CN7	568404	LEIDSCHEENDAM	NL	LSD0620	EURO
5124069	00000372661CN9	559623				EURO
5124069	00000372679CN1	553319	SINT-NIKLAAS	BELGIE	110014CT04	EURO
5124069	00000372688CN2	548725	LONDON	EG	00888CTS01	EURO
5124069	00000372825CN0	537844	HAIFA	IL	IL1086CT01	EURO
5124069	00000372882CN1	576472	MONTAUBAN		M1322104	EURO
5124069	00000372920CN9	559951	BRIIS SOUS FORGES	FRANCE	M9024408	EURO
5124069	00000372927CN4	561552				EURO
5124069	00000372951CN4	560595				EURO
5124069	00000372953CN0	564608				EURO
5124069	00000372975CN3	571486				EURO
5124069	00000374190CN7	580138				EURO
5124069	00000374531CN2	560600				EURO
5124069	00000374580CN9	540448				EURO
5124069	00000374990CN0	571091				EURO
5124069	00000375570CN9	548157				EURO
5124069	00000375774CN7	573014				EURO
5124069	00000375794CN5	570851				EURO
5124069	00000375859CN6	570722				EURO
5124069	00000375860CN4	570439				EURO
5124069	00000375867CN9	548155				EURO
5124069	00000375919CN8	561786				EURO
5124069	00000376028CN7	561644				EURO
5124069	00000376031CN1	547116				EURO
5124069	00000376067CN5	571177				EURO
5124069	00000376068CN3	547717				EURO
5124069	00000376072CN5	580933				EURO

Table 1

Purpose

The purpose of this FMI is to update the VCT installed base with the current version of software and documentation.

Furnished Materials

FMI 25372 Collector 5145026

Qty	Description	Part #
1	FMI 25372 Document	5144790-100
1	VCT M3A Interchange opt SW Key	5150379
1	Documentation Collector M3A	5136582-2
1	VCT M3A 05MW31.6 software collector	5149027
1	VCT 64 digital module config manager utility CD	5148795

Tool Required

Qty	Description	SW Version	Part #
1	FMI Document	05MW31.6	5144790-100

Personnel Requirements

"8" hours labor and "1" hours of travel

Procedure

1. Preparation or Before you Begin

- Save system state.
- Save the customer's protocols on a separate side of the state DVD.

2. LFC, VCT M3A (05MW31.6) Software Load Procedure

- Use the service CD included in the FMI kit to install VCT M3A (05MW31.6) system software.

3. VCT Field Detector Upgrade

Approximate time to complete: 15 minutes

Overview:

The detector module on board eeprom information needed for system history tracking and module identification was not programmed into the detector modules flash memory for the early systems. This procedure will update the currently installed detectors to have the correct manufacturing information.

The system software expects to see this information for module history tracking and module replacement tracking. This procedure must be run after installing the new software to maintain full compatibility.

If any modules have been replaced between system install and today, the upgrade utility will maintain the correct internal information for that module.

NOTE: After running this update procedure, do not use previously created Save State DVDs.

NOTE: For any problems during the upgrade, first exit and rerun the upgrade utility and then if problems still exist contact Eric Stepanovich for assistance 262-424-0136.

NOTE: This procedure must **NOT** be run on Application Software prior to **05MW31** software as seen from the Common Service Desktop.

Procedure:

1. Check the software revision listed on the Common Service Desktop. This procedure must not be run on Application Software prior to **05MW31** software revision.
2. Remove the gantry right side cover and disable Axial Drive.

3. Read the "SK" number from the detector rating plate found on one end of the detector carbon window as seen through the scan window gap in the covers. Pull the scan window off the gantry covers if needed.

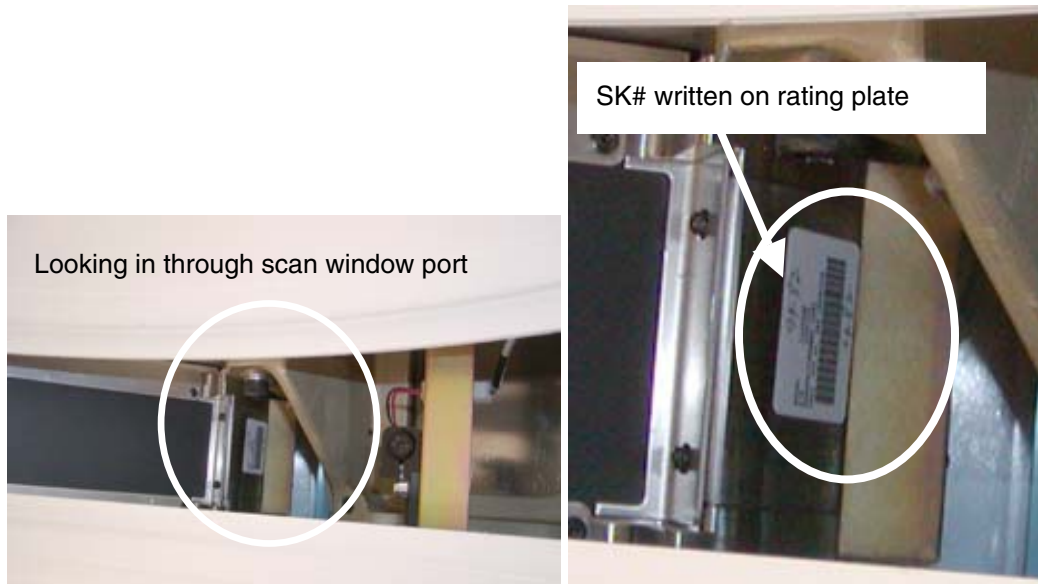


Figure 1 Detector rating plate location

4. Put the scan window back in place and enable axial drive.
5. Run DasTools maRatio's test to confirm the current list of bad pixels reported to the gesyslog.

Example gesyslog message for pixels that are mapped out:

Thu Aug 4 07:39:13 2005

Host: bay45
230007009

Proc: badchanmgr

Error:

File: BadChannelMap.cxx

Line: 742

Exclusion Test - Row 1:32A, Channel 689, Module 44 Short to ground* PASSED
CHECKS *.

6. Insert the VCT Digital Module Configuration Manager utility CDROM (5148795) into the host CD drive.
7. Open a Unix shell and type the following commands:

mount /mnt/cdrom

java -jar /mnt/cdrom/frmMenu.jar

NOTE: This must be performed as ctuser that is the default user. Do not change to root or any other user.

NOTE: The Unix shell will show what the tool is doing during the process. Watch the messages there as any problems during the upgrade can be seen in the Unix shell. Message examples would be Reading, Writing, and Validating each Module. Having the output stop with an error can usually be recovered by exiting and restarting the upgrade utility using the java command in this step.

8. The Unix shell will indicate that the system is backing up the current detector config. When the backup is done the following GUI will appear as shown in Figure 2. Select the **"Set Detector Number"** button at the top of the GUI.

The system will read the current detector and the window shown in Figure 3 will appear.

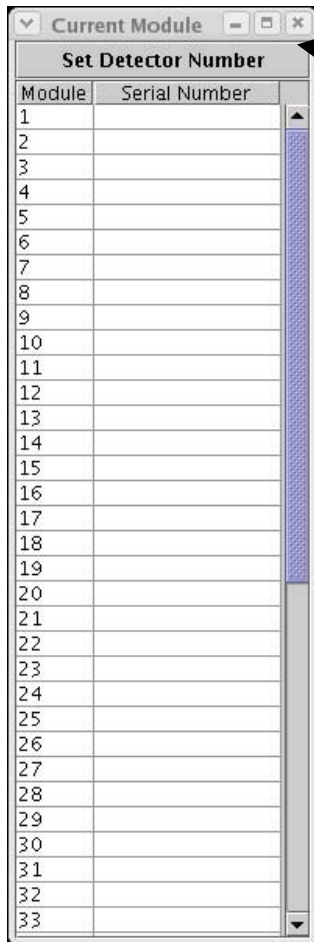


Figure 2 GUI example

"Set Detector Number" button

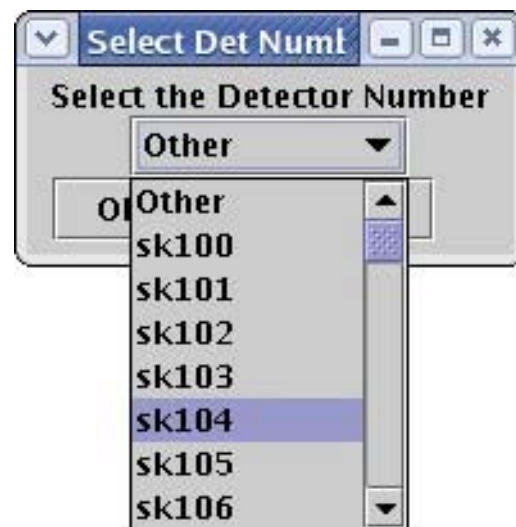


Figure 3 Detector SK# menu

9. After the modules have been read, the popup window shown in Figure 3 will appear. From the pull down menu shown in Figure 3 select the **"SK#"** that matches the value written on the detector rating plate.

NOTE: If your SK# is not on the list use the **"Other"** entry.

10. Select **"OK"** from the window and confirm your detector SK selection in the acknowledgement window that appears. The upgrade utility will now read the current detector modules and update the onboard flash memory with the correct information. The update activity can be seen in the Unix shell that was used to start the upgrade utility. When done updating the detector, the utility will automatically run a gantry hardware reset. This will trigger the system to read the modules and update the appropriate host configuration tracking files.

- a. If the update utility is interrupted in progress of updating the detector, get the system back up and running and start the update utility again.
 - b. If any java errors show up during the update process, restart the tool and try again.

Exception (Ignore this if seen): Module X: has invalid data, continuing to next a2d.
 - c. If there is a read, write, or verification error seen in the Unix shell, the utility will try a second time automatically. If the second pass fails the utility will let you know there is a problem. Cycle power to the DAS and try the upgrade again.
 - d. If the gantry hardware reset does not complete successfully after the detector modules are updated, try the hardware reset manually from the "System Resets" selection on the Common Service Desktop GUI or reset power using the gantry service switch panel. Troubleshoot as necessary. When the Gantry reset is able to complete successfully, rerun the Detector upgrade utility.
11. When the message "Detector Upgrade Complete!!" shows up in the Unix shell, the upgrade is done.
 12. To exit the GUI select the box in the upper left and select close to exit.
 13. Type the following command to unmount the CD:
umount /mnt/cdrom
 14. Remove the Upgrade CD.
 15. Run DasTools maRatio's test to confirm the list of bad pixels reported to the gesyslog.
 - a. If the list is different than what was reported to the gesyslog prior to the detector data upgrade, call Eric Stepanovich for assistance 262-424-0136.

NOTE: A bug has been seen for the gesyslog viewer right after running the detector upgrade utility. If the reported bad channels look to be a list of the same 1 or 2 channels instead of all different channels, then look you can look at the gesyslog from a Unix shell or run maRatio's again. The problem is in the gesyslog veiwer since the data in the log itself is correct. A second pass of maRatio's displays fine. After a system reboot, the log will also show the correct entries.

Finalization:

No finalization steps required.

4. Full system calibration

After the a successful software load, complete the following calibration procedures:

- Collimator Calibration
- Detailed Calibration
- FastCal

5. Save State Configuration file edit – *Must be completed before Save State is performed.*

ISSUE: Save System State does not save the necessary detector configuration files for 05MW31.6.V40_H_V64_G_GTL software release. This is a known software bug, which will be fixed in the M4 software release. Failure to perform the defined corrective action will cause FRDM Wizard Tool failures when replacing detector modules.

CORRECTIVE ACTION: Create a backup copy and edit the Save System State configuration file to correct a mistake in a filename. This will then have the Detector configuration file saved with each Save State. Make sure a Save State is performed after this update to have a complete State backup.

EFFECTIVITY: This applies ONLY to 05MW31.6.V40_H_V64_G_GTL and must be performed after each completed software Load. Keep a copy of this Save State Configuration file edit procedure with your State DVD as a reminder for any subsequent LFC.

PROCEDURE START:

Step 1:

Open a shell and execute the following commands: (Your entries/mouse actions are in bold type)

```
{ctuser@hostname} [1] cd /usr/g/config
```

```
{ctuser@hostname} [2] ls -al Sys*
```

See Figure 1 for output


```
{ctuser@hostname} [3] cp SysState.cfg.hp SysState.cfg.hp.lfc
```

```
{ctuser@hostname} [4] ls -al Sys*
```

See Figure 1 for output

```
{ctuser@hostname} [5] nedit SysState.cfg.hp
```

New editor window opens. See Figure 2.

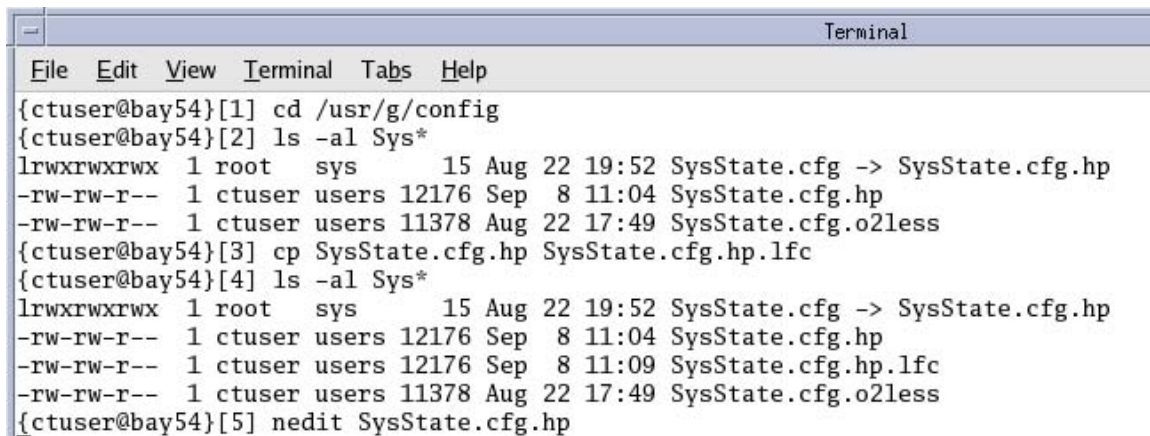

```

Terminal
File Edit View Terminal Tabs Help
{ctuser@bay54}[1] cd /usr/g/config
{ctuser@bay54}[2] ls -al Sys*
lrwxrwxrwx 1 root sys 15 Aug 22 19:52 SysState.cfg -> SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 12176 Sep 8 11:04 SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 11378 Aug 22 17:49 SysState.cfg.o2less
{ctuser@bay54}[3] cp SysState.cfg.hp SysState.cfg.hp.lfc
{ctuser@bay54}[4] ls -al Sys*
lrwxrwxrwx 1 root sys 15 Aug 22 19:52 SysState.cfg -> SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 12176 Sep 8 11:04 SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 12176 Sep 8 11:09 SysState.cfg.hp.lfc
-rw-rw-r-- 1 ctuser users 11378 Aug 22 17:49 SysState.cfg.o2less
{ctuser@bay54}[5] nedit SysState.cfg.hp


```

Figure 1

Step 2:

In the editor window **Scroll Down** until the following line is seen;

```
charactFileList_46 = { host = OC, file = /usr/g/service/state/config_tracker/masterDCF.xml, req = N}
```

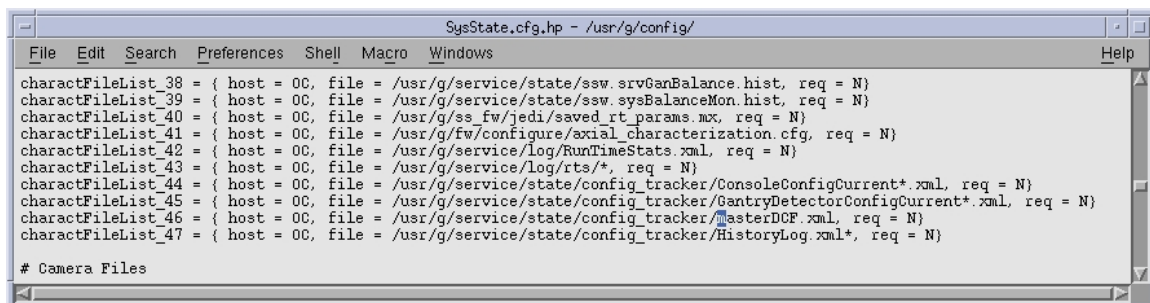

```

SysState.cfg.hp - /usr/g/config/
File Edit Search Preferences Shell Macro Windows Help
charactFileList_38 = { host = OC, file = /usr/g/service/state/ssw.srvGanBalance.hist, req = N}
charactFileList_39 = { host = OC, file = /usr/g/service/state/ssw.sysBalanceMon.hist, req = N}
charactFileList_40 = { host = OC, file = /usr/g/ss_fw/jedi/saved_rt_params.mtx, req = N}
charactFileList_41 = { host = OC, file = /usr/g/fw/configure/axial_characterization.cfg, req = N}
charactFileList_42 = { host = OC, file = /usr/g/service/log/RunTimeStats.xml, req = N}
charactFileList_43 = { host = OC, file = /usr/g/service/log/rts/*, req = N}
charactFileList_44 = { host = OC, file = /usr/g/service/state/config_tracker/ConsoleConfigCurrent*.xml, req = N}
charactFileList_45 = { host = OC, file = /usr/g/service/state/config_tracker/GantryDetectorConfigCurrent*.xml, req = N}
charactFileList_46 = { host = OC, file = /usr/g/service/state/config_tracker/masterDCF.xml, req = N}
charactFileList_47 = { host = OC, file = /usr/g/service/state/config_tracker/HistoryLog.xml*, req = N}

# Camera Files


```

Figure 2

Using the cursor, **left click** and highlight the first lower case "m" in "masterDCF.xml" only. See figure 2.

With this letter highlighted type a capital M on the keyboard: **Shift + M**

The lower case "m" has been changed to an uppercase "M". See figure 3.

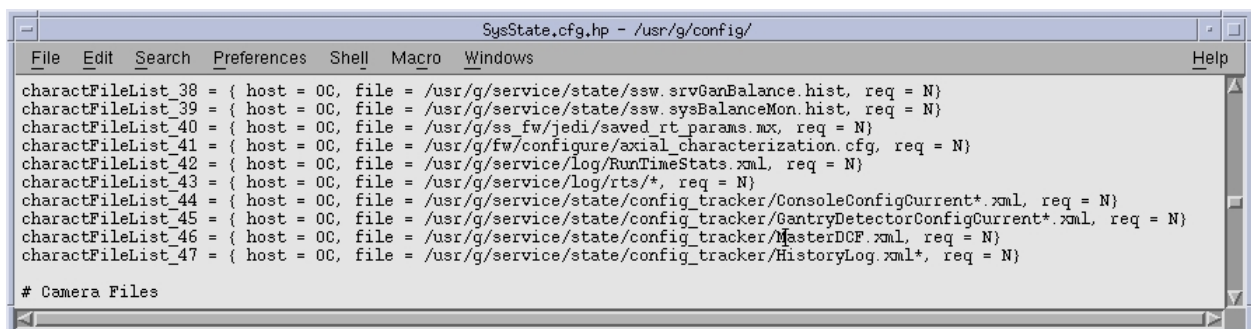

```

SysState.cfg.hp - /usr/g/config/
File Edit Search Preferences Shell Macro Windows Help
charactFileList_38 = { host = OC, file = /usr/g/service/state/ssw.srvGanBalance.hist, req = N}
charactFileList_39 = { host = OC, file = /usr/g/service/state/ssw.sysBalanceMon.hist, req = N}
charactFileList_40 = { host = OC, file = /usr/g/ss_fw/jedi/saved_rt_params.mtx, req = N}
charactFileList_41 = { host = OC, file = /usr/g/fw/configure/axial_characterization.cfg, req = N}
charactFileList_42 = { host = OC, file = /usr/g/service/log/RunTimeStats.xml, req = N}
charactFileList_43 = { host = OC, file = /usr/g/service/log/rts/*, req = N}
charactFileList_44 = { host = OC, file = /usr/g/service/state/config_tracker/ConsoleConfigCurrent*.xml, req = N}
charactFileList_45 = { host = OC, file = /usr/g/service/state/config_tracker/GantryDetectorConfigCurrent*.xml, req = N}
charactFileList_46 = { host = OC, file = /usr/g/service/state/config_tracker/MasterDCF.xml, req = N}
charactFileList_47 = { host = OC, file = /usr/g/service/state/config_tracker/HistoryLog.xml*, req = N}

# Camera Files

```

Figure 3

Save the file and exit:

Left Click the **File** Pull down and select **Save**.


Left Click the **File** Pull down and select **Exit**.

Step 3:

In the open shell window type the following;

```
{ctuser@hostname} [6] diff SysState.cfg.hp SysState.cfg.hp.lfc
```

Verify the only difference between these two files is the "M" in MasterDCF.xml is capitalized.
See Figure 4.

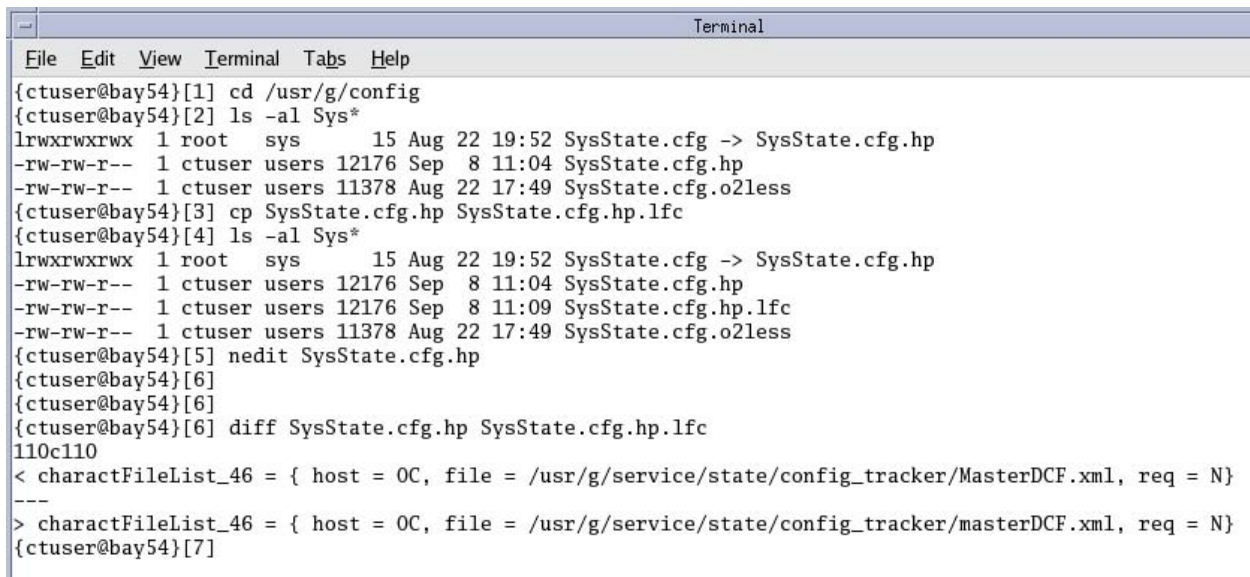

```

File Edit View Terminal Tabs Help
{ctuser@bay54}[1] cd /usr/g/config
{ctuser@bay54}[2] ls -al Sys*
lrwxrwxrwx 1 root sys 15 Aug 22 19:52 SysState.cfg -> SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 12176 Sep 8 11:04 SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 11378 Aug 22 17:49 SysState.cfg.o2less
{ctuser@bay54}[3] cp SysState.cfg.hp SysState.cfg.hp.lfc
{ctuser@bay54}[4] ls -al Sys*
lrwxrwxrwx 1 root sys 15 Aug 22 19:52 SysState.cfg -> SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 12176 Sep 8 11:04 SysState.cfg.hp
-rw-rw-r-- 1 ctuser users 12176 Sep 8 11:09 SysState.cfg.hp.lfc
-rw-rw-r-- 1 ctuser users 11378 Aug 22 17:49 SysState.cfg.o2less
{ctuser@bay54}[5] nedit SysState.cfg.hp
{ctuser@bay54}[6]
{ctuser@bay54}[6]
{ctuser@bay54}[6] diff SysState.cfg.hp SysState.cfg.hp.lfc
110c110
< charactFileList_46 = { host = 0C, file = /usr/g/service/state/config_tracker/MasterDCF.xml, req = N}
---
> charactFileList_46 = { host = 0C, file = /usr/g/service/state/config_tracker/masterDCF.xml, req = N}
{ctuser@bay54}[7]

```

Figure 4

Step 4:

Perform a Save system state "All".


6. Collimator ISO adjuster screws

Problem

Screws attaching the ISO adjuster block to the collimator interposer plate may be missing Loctite, which could cause the screws to loosen or fall out.

Solution

- Remove the screws attaching the adjuster block to the interposer plate, one screw at a time.
- Apply Loctite 242 (blue) and reinstall screws. Torque to 2.9 NM. See figure.

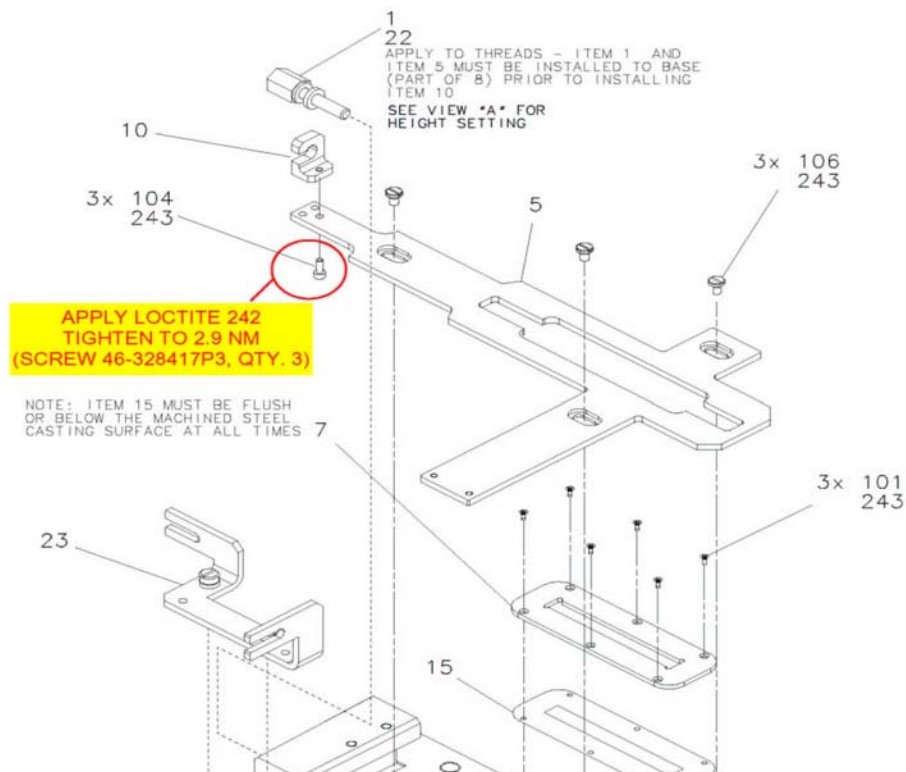

7. TGPU Board Verification

- Verify TGPU board Revision

Look at the board part or use the following command:

```
[root@ct03]# rsh tgp
```

```
->print_cpu_board_revision
```

PartNo: **2349697-03C4 (Variation:3** Revision:4868)

Value = 49 = 0x31 = '1'

```
->logout
```

If your TGPU is not a –03, "xxxxxxx-03xx (**Variation:x** Revision:xxxx)" you will need to e-mail:

Eric.Stepanovich@med.ge.com

There will only be a very small percentage of circuit boards that are not Rev – 03.

8. WW/WL Presets

After Restore System State and running the cleanup script (postRestoreCleanup.vctM3A).

Go to Image Works

Start Viewer

Click "User Preferences"

Click "Save As Default"

Wait about 10 - 15 seconds to allow the update to complete.

Click "Browser"