

**ADMINISTRATIVE INFORMATION****FMI No: 25374****FMI Title: VCT Failsafe Collimator Upgrade**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS - AM:	Issue Date: Feb. 3 rd , 2006
☒ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date: Jun. 30 th , 2006

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.:** 5150103-100**Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas

GEMS - Europe

G.E. Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe
FMI Admin. 322
BP34
F78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25374

FMI Title: **VCT Failsafe Collimator Upgrade**

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments _____

Effectivity

Existing VCT system shipped prior to Q2'05.

Univ. of Washington	206598VCT2
Charite	CC0618CT07
AsiaMedic	0847210010
Toki Municipal Hospital	QS0001
Handai Fuzoku (Osaka)	QX0008
Dortmund	HC0074CT04
Huntsville Medical Center, AL	256MAINVCT
Cardiology Associates, AL	251340VCT
West Suburban, IL	630972VCT
CHUV Lausanne	260010CT15
RIKSHOSPITALET HF	N01016CT06
University of Michigan	734936VCT
UNIVERSITA' DEGLI STUDI DI ROMA TOR VERGATA, Italy	A5333029
BEIJING PLA NO.301 HOSPITAL	082421040148
AMI	609652VCT
Beijing Xuanwu Hospital	082421040151
Miami Baptist, FL	305596VCT
Boston Med Ctr, MA	617414VCT1
North Shore Univ., NY	516562VCT
Bronson Methodist Hospital, MI	269341VCT
Mobile Infirmary	251435CT3
IMorristown, NJ / Mountainside	973429VCT
East River Medical imaging, NY	212288CT64
St. Mary's Hospital, MI	989497CT
Methodist Hospital, TX	713441MCT9

Purpose

The purpose of this FMI is to inspect the VCT system and determine if the existing collimator is of FAILSAFE type or not. If not, the existing collimator must be replaced with the FAILSAFE version included in the FMI kit.

NOTE: All VCT systems shipped prior to Q2'05 should have the NON-FAILSAFE version.

Prerequisite FMI / Related FMIs

None

Furnished Materials

Qty	Description	Part #
1	FMI Document	5150103-100
1	Failsafe Collimator	5130001

Tools Required

Qty	Description
1	Standard FE service kit.


Personnel Requirements

"8" hours labor and "1" hours of travel.

Procedure

Preparation

1. Verify that the existing collimator in the system is the FAILSAFE version by looking through the Mylar window. Part number will be (5130001).

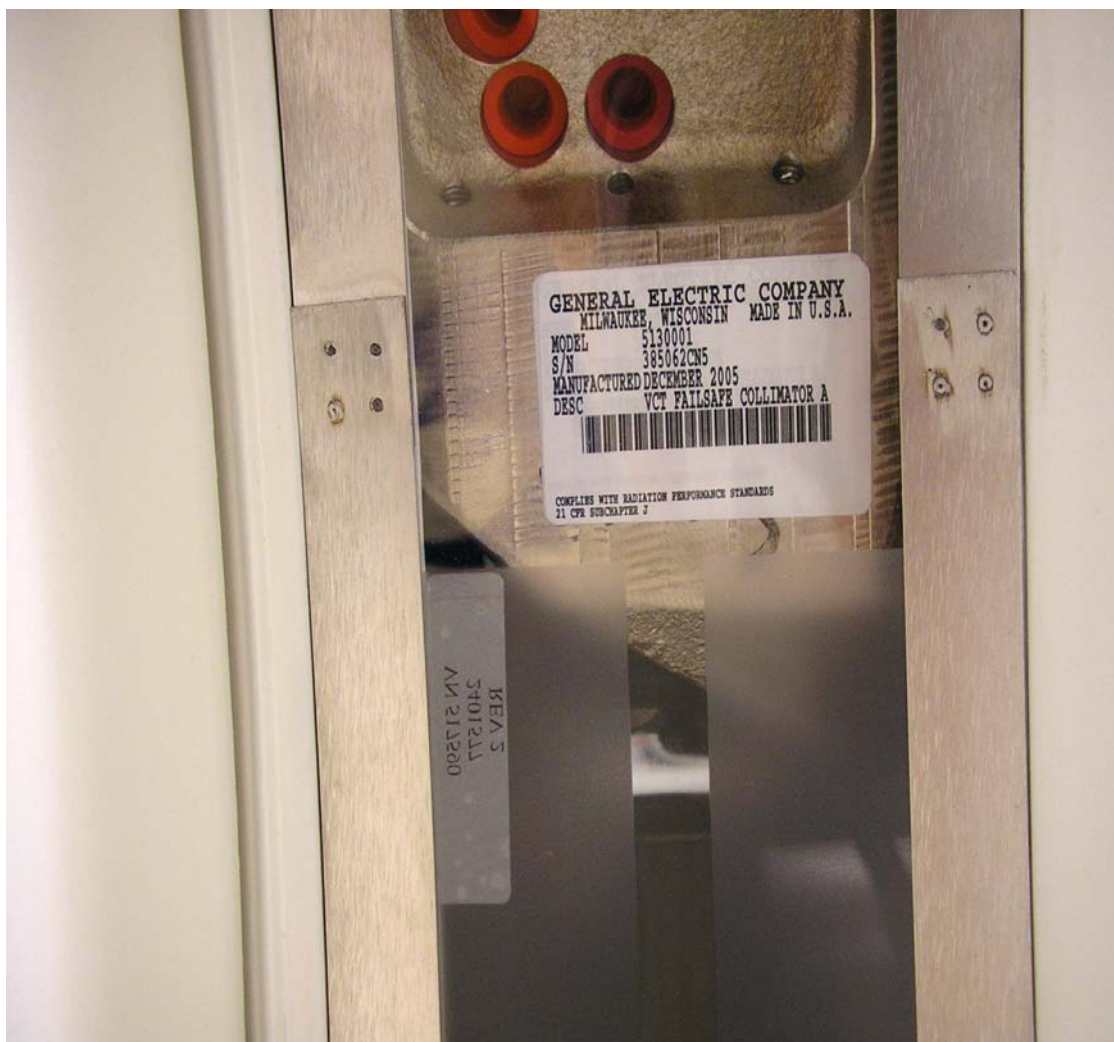

View through the Mylar window.



View when looking from right side, tube @ 2:00.

2. If the existing collimator in the system is NOT FAILSAFE replace the collimator with the FAILSAFE version included in the FMI kit. Go to the "**Collimator Replacement**" section on page 5.
3. If the existing collimator in the system is FAILSAFE stop here. The FMI does not need to be performed. Return unused FMI kit to:

HM Product Solutions, LTD.

3225 Gateway Road, Suite 250

Brookfield, WI 53045

262-754-0424

Collimator Replacement

1. Please refer to service CD **5116423 - LS 7.X Service Methods, Replacement tab - Gantry - collimator assembly replacement**, follow accordingly. The tube replacement procedure is linked inside the collimator replacement procedure.

Finalization

1. From the **Common Service Desktop, Calibration Menu** select [**Gantry Service balance**] and perform a balance check. Rebalance the gantry if necessary. Should only require balancing if gantry was almost out of specification to start.
2. Perform a test scan using the 20cm QA phantom to verify system operation.
3. Close the FMI.