

**ADMINISTRATIVE INFORMATION****FMI No: 25375****FMI Title: CT H1 Gantry Helios Collimator Filler Bar Front Cover Interference Kit**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

- ☒ **Staged FMI** **Customer List** provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
- ☒ **Parts** and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
- ☐ **Parts Ordering:** GEMS - Americas, Contact your Region logistics coordinator
GEMS - Europe, Order FMI kits from GPO/EPL Buc
GEMS - Asia, Contact your region logistics coordinator
- ☐ **No Parts** are required to complete this FMI.
- ☒ **Batch with** FMI 25360
- To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
- ☐ **Non-Staged** **No Customer List is provided.** Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked: 2180482****FMI Instructions Part No.: 5159481-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25375

FMI Title: CT H1 Gantry Helios Collimator Filler Bar Cover
Interference Kit

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments _____

Effectivity

All CT H1 gantries, tracked gantry model number 2180482.

Purpose

This FMI kit provides the parts and instruction to eliminate a mechanical interference condition between the CT H1 gantry front cover laser alignment windows and the Helios collimator filler bar deployed by FMI 25360 for 2214768-3 and 2214768-4 collimators.

The mechanical interference problem is resolved by replacing the laser windows and brackets with laser windows that are glued directly to the front gantry cover.

Prerequisite FMIs

None

Related FMIs

FMI 25360, Helios Collimator Xray Scatter Filler Bar for 2214768-3 & -4 Collimators (p/n 5133594)

Furnished Materials

FMI 25375 Kit Part Number = 5159480			
Qty	Description	Part #	Comments
1	FMI Document	5159481-100	
1	Adhesive cartridge, 3M DP100 Structural	5130489	3-5 minute work time. Hardening time increases with amount applied (estimate = 7 minutes).
1	Gun, Applicator for DP100 Adhesive	5147976	
6	Nozzle, Mixing tip for applicator gun	5147983	
2	Window, laser aperture	5150546	

Tools Required

No additional tools; standard FE tool kit required

Rubber gloves from CT system installation kit at site


Personnel Requirements

1 FE	resource
1.5 hour	labor
0 hour	travel (batch with FMI 25360, Helios Collimator Xray Scatter Filler Bar)

Procedure

1. Position table and gantry to gain the greatest access to the front cover laser windows located at 3 and 9 o'clock positions.
Open front cover as far as possible.
2. Turn off gantry power and perform LOTO steps for personal safety.
3. Manually rotate the x-ray tube/collimator to the 5 o'clock position to prevent dirt from falling onto it.
(View from front of scanner. The gantry/collimator may be rotated by hand – manually).
4. Cover the collimator to prevent dirt from falling onto it.
Use the work pad found in your service tool kit as a cover.
5. Mark the location of both laser windows on the inside of the cover to help locate the “new” laser window provided in this kit.


6. Remove the two laser window / bracket assemblies. Refer to Picture 1.
Remove screws and window/bracket, and discard.

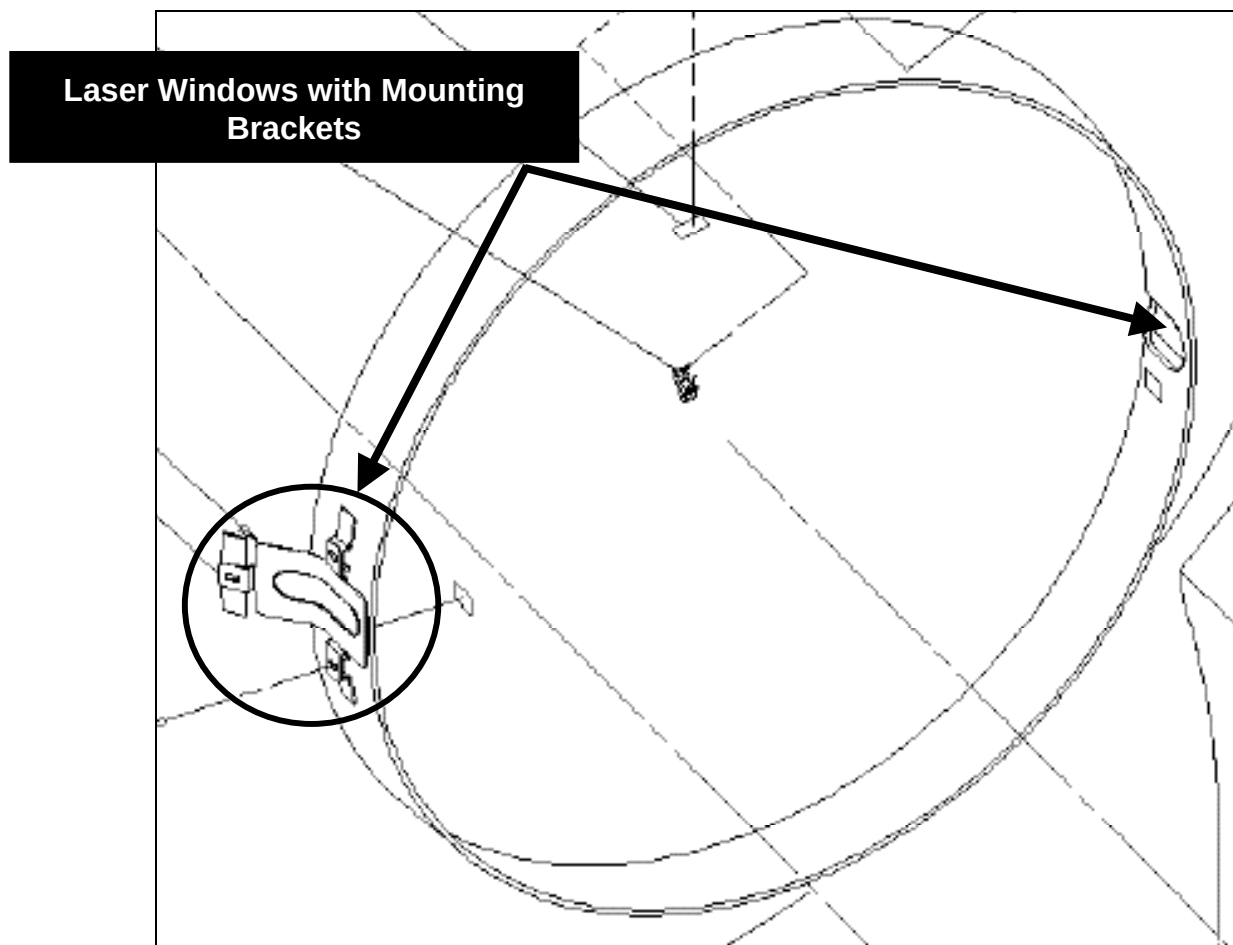

Picture 1, Remove Laser Windows and Brackets from Front Gantry Cover


7. Remove the mounting tabs from the front gantry cover.
Refer to Picture 2 for an illustration of mounting tabs to be removed.
- Use pliers to grab / twist / remove the tabs from the fiberglass.
Three tabs exist per bracket. Only remove two tabs, located at 6 & 12 o'clock positions.
 - Remove the fiberglass material that held the mounting tabs.
Use needle nose pliers to grab/twist/remove the fiberglass (or flat blade screwdriver).

NOTE: Removing these mounting tabs creates a mechanical clearance to prevent the collimator filler bar installed by FMI 25360 (Helios Collimator Xray Scatter Filler Bar) from interfering with the front cover on H1 gantries.

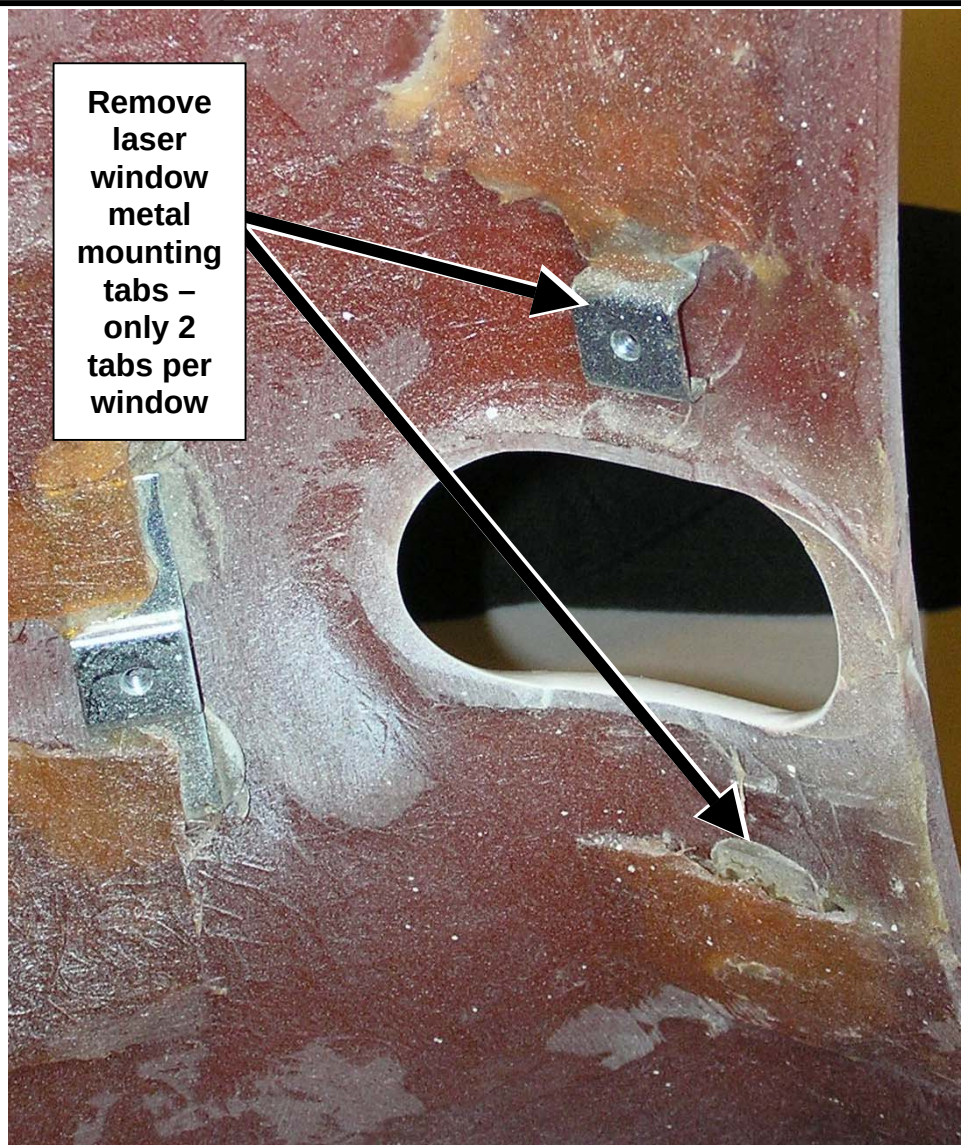

Picture 2, Mounting Tabs to Remove – Qty 2 per Laser Window

8. If the front gantry cover is RIM material, then roughen the inside, painted surface to remove paint to create a better adhesion condition for gluing the new laser window in place.
Wipe away any dust.
9. Test the position of a new laser window, centering it in the cover's laser opening.
The notched end of the window must be positioned towards the CT gantry.
Refer to Picture 3 for "notch" on window.
NOTE: Do not remove the protective tape from the laser window.

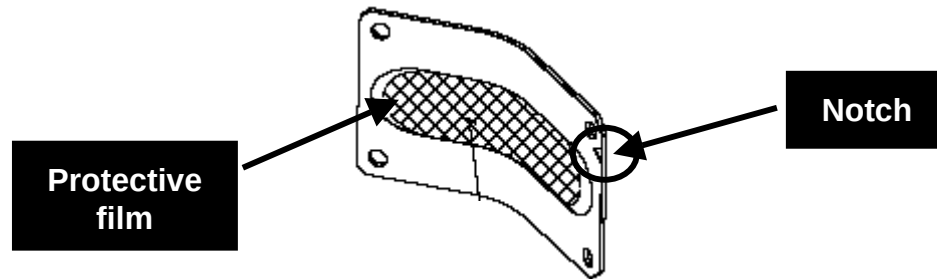

Picture 3, Notch on Laser Window

10. Make any necessary positioning marks that will help locate the window before gluing (glue work time = 3-5 minutes).

NOTE: Position the laser window with notch located towards scanner / gantry bore.

11. Setup the epoxy / glue gun as follows:

- a. Locate gun from FMI kit
- b. Locate the epoxy cartridge from the FMI kit
- c. Load the gun with the epoxy cartridge
- d. Locate the epoxy nozzle
- e. Attach the nozzle to the end of the cartridge

12. With a laser window in hand, form the laser window to fit the shape of your front cover (inside front cover)

- a. Position laser window with notch located towards scanner / gantry bore. Refer to Picture 3 for view of laser window notch.
- b. If needed, flex / shape the window to conform it to the shape of your cover (inside cover).
- c. Confirm the fit, re-position the laser window to the inside of the front cover

13. Apply a glue bead to the gantry cover's laser opening, inside surface, corners only.

14. Apply a glue bead to the roughened surface of the laser window.

Do not apply outside the roughened area.

15. Position the laser window to the inside of the cover, using the positioning marks.

16. Hold (secure) the laser window in place until the glue hardens.

NOTE: Work time = 3 - 5 minutes. Time to harden increases with amount applied.

7 minutes is approximately the time to harden.

17. Replace the nozzle on glue cartridge before installing the next laser window.

NOTE: Release plunger pressure on glue cartridge before replacing nozzle.

18. Repeat steps to position, glue, and hold the other laser window (located at 9 o'clock position).

19. Again, replace the nozzle on glue cartridge before filling the gap between the laser window and gantry cover with glue to create a fluid-tight seam.

NOTE: Release plunger pressure on glue cartridge before replacing nozzle.

20. Return to the laser window located at the 3 o'clock position.

Fill the gantry cover to laser window gap with glue to prevent fluids from entering the system.

Fill gap at bottom of the laser window by filling the gap from the front of the cover.

Fill gap at top of the laser window by filling the gap from the inside of the cover.

Fill the four laser window corner holes to improve the bond to the front cover.

NOTE: Do not apply an excessive amount of epoxy at one time – it may run.

Work between the two laser windows when filling the gap. This will allow the small bead of glue to harden between applications, preventing glue runs and glue hardening in the tip.

21. Repeat the gap-filling steps for the laser window located at the 9 o'clock position.

22. If the gantry cover is fiberglass, apply glue to the fiberglass areas where the mounting tabs were removed to encapsulate any exposed fibers.

Use the applicator tip to spread the glue across this area.

Do not use your finger to spread the glue on the fiberglass surface.

23. With clean hands, remove the protective film from both sides of the laser window located at 3 o'clock..

24. With clean hands, remove the protective film from both sides of the laser window located at 9 o'clock.

25. Proceed to FMI 25360 to install the collimator filler bar per FMI 25360's instruction.

Finishing Up

Perform the necessary steps to close dispatches for FMI 25375 and 25360.