

**ADMINISTRATIVE INFORMATION****FMI No: 25385****FMI Title: CT Gantry Control Touch Panel Replacement**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS - AM:	Issue Date:
☒ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

- ☒ **Staged FMI** **Customer List** provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
- ☒ **Parts** and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
- ☐ **Parts Ordering:** GEMS - Americas, Contact your Region logistics coordinator
GEMS - Europe, Order FMI kits from GPO/EPL Buc
GEMS - Asia, Contact your region logistics coordinator
- ☐ **No Parts** are required to complete this FMI.
- ☐ **Batch with** _____
To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).
- ☐ **Non-Staged** **No Customer List is provided.** Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS – Americas	GEMS - Europe	GEMS – Asia
GE Medical Systems Insite Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:

The issue being corrected by this FMI has been addressed in Spare Part Stock. When performing this FMI, one needs to check if the Gantry Control Panels have already been replaced with new version. See instructions for parts numbers.

If Gantry Control Touch Panels are not needed, return the FMI kit to HM Product Solutions Ltd. and code FMI Service Dispatch as Code 2 (FMI Not Required). See instruction for shipping information.

Unit Tracked: Refer to Effectivity list on page 3.**FMI Instructions Part No.: 5198284-100, Rev 2****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas

GEMS - Europe

G.E. Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe
FMI Admin. 322
BP34
F78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25385

FMI Title: CT Gantry Control Panel Replacement

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments

Table of Contents

Table of Contents	3
Revision History	3
Effectivity	4
Purpose	5
CUSTOMER IMPACT	5
Prerequisite FMIs	6
Related FMIs	6
Furnished Materials	6
Tool Requirements	6
Personnel Requirements	6
Before You Begin	7
Safety	7
Gantry Control Touch Panel Replacement	7
Procedure	8
GANTRY COVER REMOVAL	8
REMOVING GANTRY CONTROL TOUCH PANEL	9
INSTALLING NEW CONTROL TOUCH PANEL	10
Finalization	11

Revision History

Revision	Date	Reason for Change
1	4/30/2007	Initial release. Updated from HII Review. Initial Engineering review for milestone - F2
2	6/14/2007	Pilot Feedback

Effectivity

The effectivity for this FMI applies to all of the LightSpeed 5.X-7.X systems manufactured between April 1, 2004 and March 31, 2006. The models tracked in effectivity are below:

Product	Gantry Assembly
VCT/DDAS64	5124069
VCT/DDAS32	5126093
HP100/MDAS16	2357739
HP100/GDAS16	2357739-4
HP100/GDAS16/SKD	2357739-6
HP100/GDAS16/EMC2	2357739-7
HP100/GDAS16/EMC2/ASTEC	2357739-8
RT/MDAS8	2374681
RT/MDAS4	2374681-2
RT/GDAS4/EMC2/ASTEC	2374681-4
HP60/GDAS16	2378661-3
HP60/GDAS16/EMC2	2378661-10
HP60/GDAS16/EMC2/ASTEC	2378661-11
HP60/GDAS8	2383549-2
HP60/GDAS8/EMC2	2383549-10
HP60/GDAS8/EMC2/ASTEC	2383549-11
HP60/GDAS4	2383550-2
HP60/GDAS4/EMC2	2383550-10
HP60/GDAS4/EMC2/ASTEC	2383550-11

Note: PET DLS or DST systems are not included in this FMI Effectivity.

PET Gantry Control Touch Panels have different button assignments and labeling.

Purpose


Customer Impact

1. iTrak Reports

13041023, 13041061, 13052977, 13064537, 13064805, 13074236, 13076712,
13080311, 13082794, 13082531, 13092221, 13093123, 13093542

2. Problem: Control Panel Buttons Stop Functioning and Fall Out

- a. Type of material used on elastomer holding keypads was too brittle and fails between 100-200K cycles. Depending on use of the CT system this could be as early as 1 year.
- b. Stainless Steel dome switch (located under buttons, contact to circuit board) also lost contact effectiveness between 100-200K cycles (switch activations).

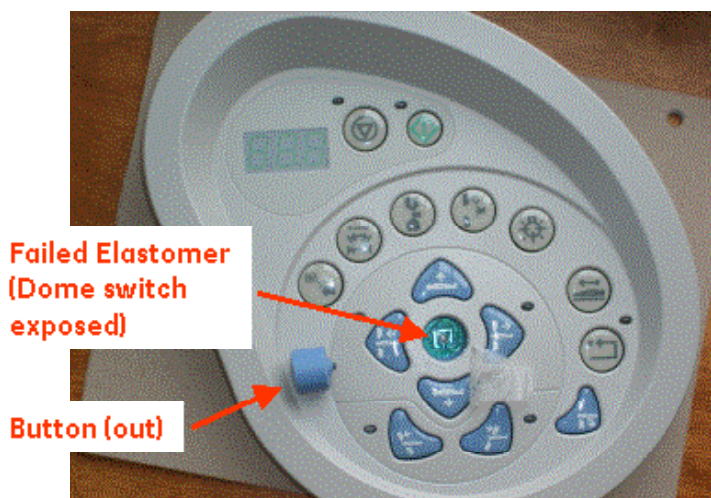

Prerequisite FMIs

None.

Related FMIs

None.

Furnished Materials

FMI Kit 5198284 consists of the following materials:

Item	Part No	Qty	Description
1	5198284-100	1	FMI 25385 Document
2	5139591	1	CT Right Single Contact, Gantry Control Touch Panel
3	5139590	1	CT Left Single Contact, Gantry Control Touch Panel

Tool Requirements

Quantity	Description
1	FE Toolkit
1	ESD wrist-band

Personnel Requirements

Resource	Labor	Travel
1 FE	1.0 hrs	1.0 hrs

Before You Begin

Before beginning this procedure, please read the gantry safety information.

Safety



⚠ WARNING

POTENTIAL FOR SHOCK.

VOLTAGE MAY BE PRESENT.

USE PROPER LOCKOUT/TAGOUT PROCEDURES BEFORE REPLACING COMPONENTS.

Gantry Control Touch Panel Replacement

This FMI supplies only two (2) replacement panels, one left and one right. The system actually has four (4) control panels, one left and one right and on both front and back of gantry. The person performing this FMI will need to decide which panels to replace. It is recommended that the most used/accessed panels be replaced, typically the two panels on the front of the gantry.

Note:

The new/improved panels have been available from spare's stock since March 2006. Chances are the system may already contain the new panels if problematic panels failed in the installed base after March 2006. Please reference the table below for the part numbers of these panels.

Return unused FMI kit to:

HM Product Solutions, LTD.

3225 Gateway Road, Suite 250


Brookfield, WI 53045 Tel: 262-754-0424

Gantry Control Touch Panel Versions	Left Side	Right Side
Problematic Gantry Control Touch Panels	2330695	2330703
New Gantry Control Touch Panels	5139590	5139591

Procedure

Gantry Cover Removal

1. Remove both Gantry side covers.

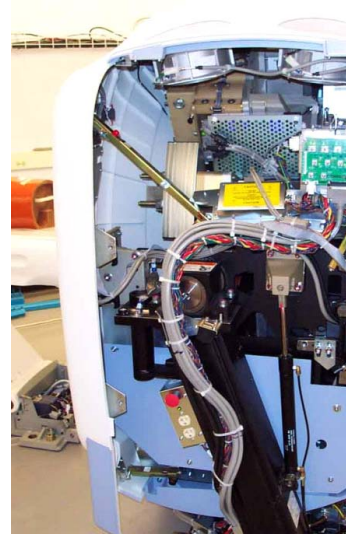

NOTICE

Always turn OFF the HVDC before the 120 VAC. Turning OFF 120 VAC power before HVDC power can result in equipment damage.



⚠ DANGER

ELECTROCUTION HAZARD.


HIGH VOLTAGE PRESENT.

USE PROPER LOCKOUT/TAGOUT PROCEDURES BEFORE REPLACING COMPONENTS.


2. Turn OFF Axial Enable, HVDC, and 120VAC Service switches to disable Gantry motion.
3. Apply Rotation Lock assembly and Lockout/Tagout procedures.

Removing Gantry Control Touch Panel

1. Remove the Gantry Control Panel. (Same for either left or right side, front or back)
 - a. Retract front cover support bracket located by the panel being replaced.

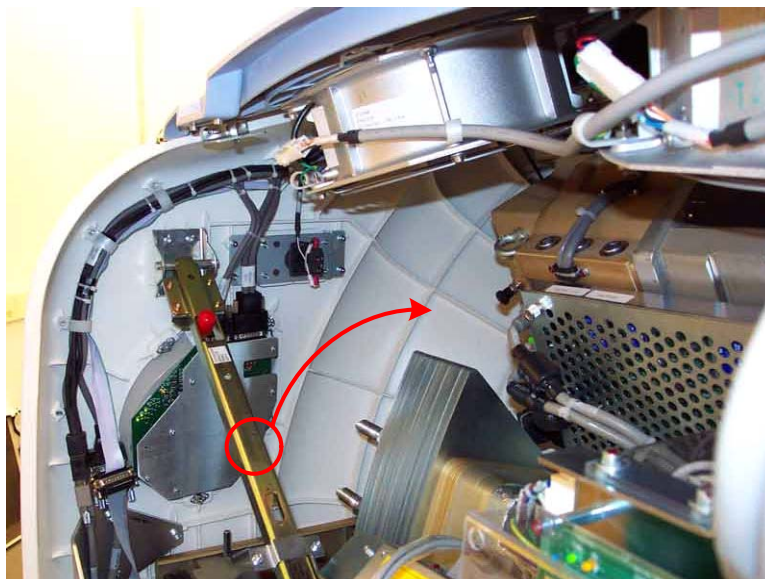

- b. Disconnect the two (2) cable D-Shell connectors from the panel.

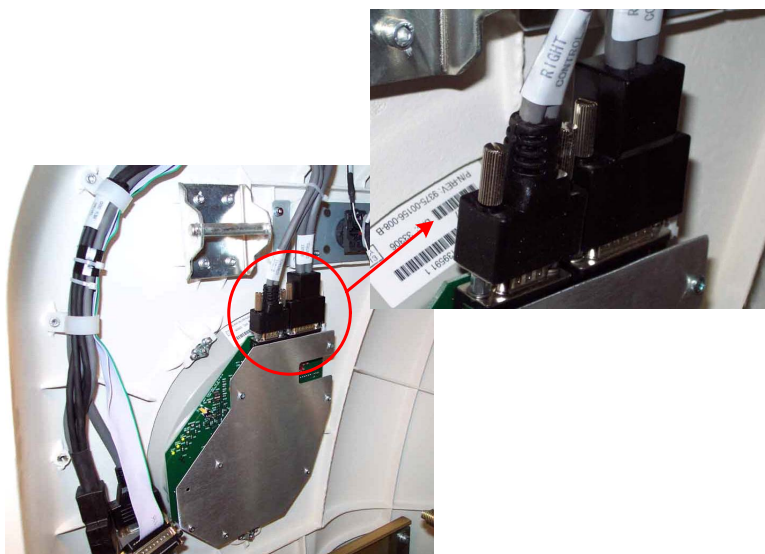

- c. Press on each ball stud till the panel releases out the Front Panel. Keep one hand on the outside of Control Touch Panel at all times to keep it from falling.



Note: The metal plate that covers the component side of the Control Touch Panel does not need to be removed. It has been removed in the picture for clarity.

2. Repeat Step 1 for the second Gantry Control Touch Panel.

Installing New Control Touch Panel

1. Install the replacement Gantry Control Panel.
 - d. Attach the Gantry Control Panel to the front cover by inserting the ball studs into the receivers.
 - e. Connect the two (2) cables to the panel.
 - f. Re-attach the front cover support bracket located by the panel being replaced.
2. Repeat Step 1 for second Gantry Control Touch Panel.

Finalization

1. Remove Rotational Lock Assembly and Lockout/Tagout procedures applied at beginning.



NOTICE

Always turn ON the 120 VAC before the HVDC. Turning ON HVDC power before 120 VAC power can result in equipment damage.

2. Turn ON 120VAC, HVDC, and Axial Enable, Service switches to enable Gantry motion.
3. Test Gantry Control Panel functionality on both sides. Test all buttons for proper operation.
4. Install the Gantry side covers.
5. Perform test scan to confirm system performance before returning system for patient scanning.
6. Close FMI service dispatch per local procedures and send old panels to Recycling for disposal.