

**ADMINISTRATIVE INFORMATION****FMI No: 25878****FMI Title: HiSpeed QX/i Ver. 100.3 S/W Upgrade**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS - AM:	Issue Date:
☒ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☐ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☐	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support		
GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2349557-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
G.E. Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201	GEMS - Europe FMI Admin. 322 BP34 F78533 Buc Cedex, France	Send to your Country FMI or Service Administrator

FMI No.: 25878

FMI Title: HiSpeed QX/i Ver. 100.3 S/W Upgrade

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer

Comments _____

Effectivity

HiSpeed QX/i systems with Gantry with the following serial numbers:

882876YM6	882869YM1	882870YM9	882871YM7	882872YM5	882868YM3
882874YM1	882878YM2	882875YM8	882873YM3	882879YM0	882877YM4

Purpose

The purpose of this FMI is to correct application software bugs on Reli Monitoring, Variviewer, and Recon (timeout).

Prerequisite FMI's

None.

Related FMI's

None.

Furnished Materials

FMI 25878 Kit (Part Number 2351135) includes the following items:

Qty	Description	Part #
1	FMI 25878 Document Set	2352212
1	Customer Letter	2350108-100
1	HiSpeed QX/i Serialized Software Set	2348418

The following items are included in the FMI 25878 Document Set

Qty	Description	Part #
1	FMI 25878 Document	2349557-100
1	LS2002 Software Installation Procedures	2348269-100

The following items are included in the Serialized Software Set

Qty	Description	Part #
1	Operating System Disk 1	2305087-2
1	Operating System Disk 2	2305088
1	Operating System Disk 3	2305089
1	Applications Software	2348417
1	LS2002 M3 Patch CD	2349518


Tool Required

Qty	Description	Part #
1	HiSpeed QX/i Class C	2344422
1	HiSpeed QX/i Class M	2344423

People Power

4 hours labor and 1 hour of travel

Steering Guide

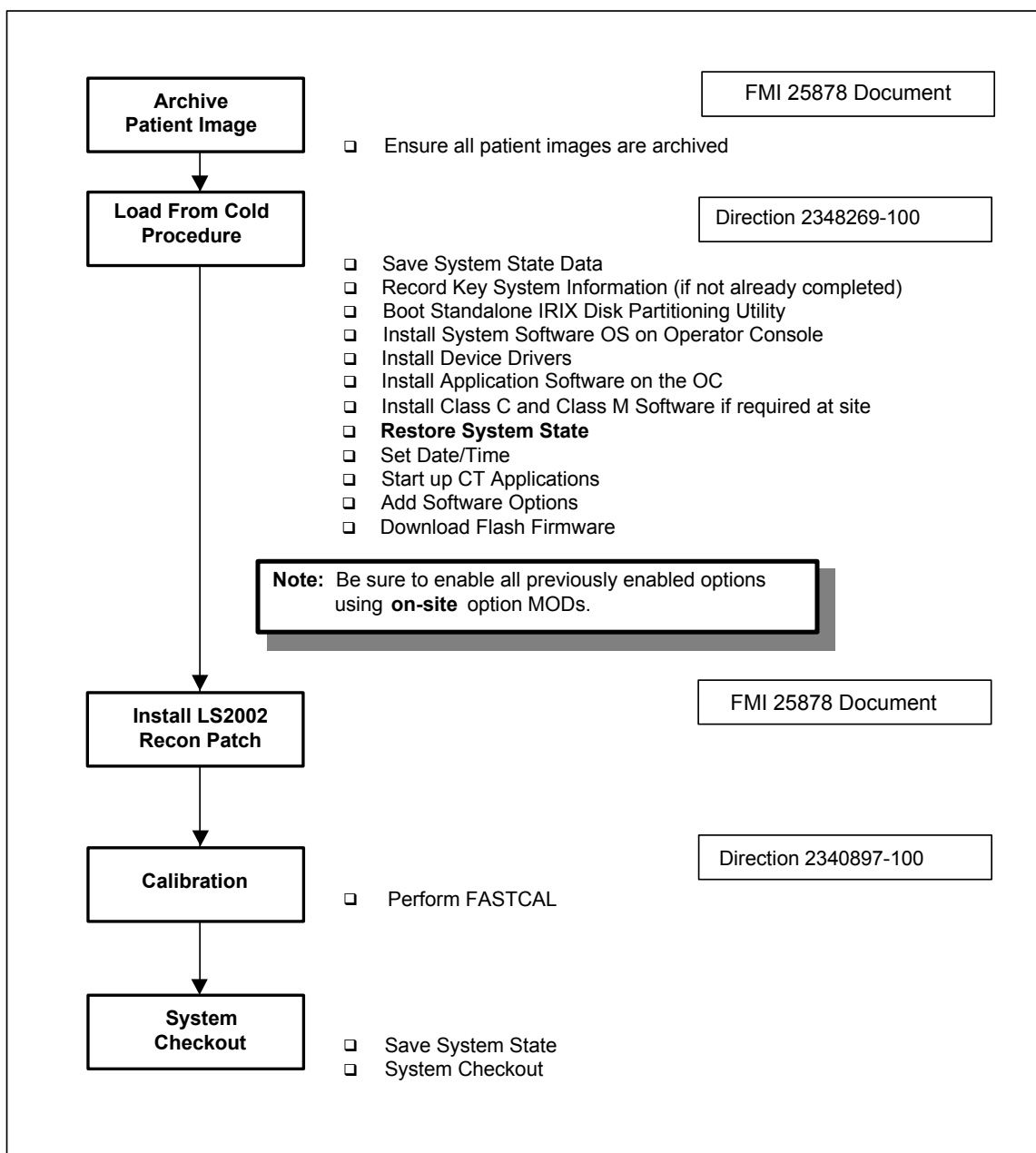

Procedure

Section 1 - Archive Patient Images

Verify that the customer has archived all patient data.

Section 2 - Software Load From Cold

Perform the Load From Cold procedure as detailed in Chapter 1 of "LS2002 Software Installation Procedures" (2348269-100).

NOTICE: It is required that a "Restore STATE" be performed following the load of Class C and Class M software. To facilitate this, it is recommended that the Class C and Class M be loaded following the Applications Load (Section 3.6).

Section 3 – Installing the LS2002 M3 Recon Patch

1. Installation:

Shut down applications by:

Service -> Utilities -> Application Shutdown

Open a UNIX shell and become root by:

type: **su -**

Enter the root password.

Insert the LS2002 M3 Recon Patch CD in the CDROM Drive:

type: **mountCD /CDROM**

type: **/CDROM/install_ls2002_m3_patch**

The install messages will appear:

.....
.....

When the installation is complete,

type: **unmountCD**

Reboot the system

2. Verification of the software patch installation:

To verify Patch is installed in OC:

type: **showprods | grep ls2002_patch_b**

The output will look similar to the following:

ls2002_patch_b	07/10/2002	Patch for ls2002 M3
ls2002_patch_b.img	07/10/2002	Patch files for recon
ls2002_patch_b .img.ocroot	07/10/2002	Files that go onto ocroot

Section 4 - Calibration

Perform FASTCAL by selecting Daily Prep on the bottom of the left monitor, then selecting FASTCAL from the displayed selections.

NOTICE: You may get a pop-up indicating: "System detects an alignment or change was performed on the System". If this happens, press the Ok button and continue with the FASTCAL.

Section 5 - System Checkout and Save State

1. Refer to the LFC Manual, and perform a Save System State. This will save the new calibrations as well as the customer and service protocols.
2. Refer to the LFC Manual, and perform the system performance checks as described in the manual. Ensure all networking and filming processes are still functional.
3. Complete and return Red Installation Product Locator Cards.

Appendix A – Uninstalling the LS2002 M3 Recon Patch

Shut down applications by:

Service -> Utilities -> Application Shutdown

Open a UNIX shell and become root by:

type: **su -**

Enter the root password.

Insert the LS2002 M3 Recon Patch CD in the CDROM Drive:

type: **mountCD /CDROM**

type: **/CDROM/uninstall_ls2002_m3_patch**

The uninstall messages will appear:

.....
.....

When the installation is complete,

type: **cd**

type: **unmountCD**

Reboot the system