



GE Medical Systems

Technical Publication

**Direction 2395230-100
Revision 0**

**5.X LightSpeed RT
Application Tips and Work-Arounds**

Copyrighted © by GE Medical Systems, 2003



Page Left Intentionally Blank

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent write "Damage In Shipment" on ALL copies of the freight or express bill BEFORE delivery is accepted or "signed for" by a GE representative or hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14-day period.

Call Traffic and Transportation, Milwaukee, WI (262) 785 5052 or 8*323 5052 immediately after damage is found. At this time be ready to supply name of carrier, delivery date, consignee name, freight or express bill number, item damaged and extent of damage. Complete instructions regarding claim procedure are found in Section S of the Policy And Procedures Bulletins. 14 July 1993.

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

IMPORTANT...X-RAY PROTECTION

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Medical Systems Group, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protection, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Medical Systems Group, its agents, and representatives have no responsibility for injury or damage, which may result from improper use of the equipment. Various protective materials and devices are available. It is urged that such materials or devices be used.

OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives. GE personnel, please use the GEMS CQA Process to report all omissions, errors, and defects in this publication.

Page Intentionally Left Blank

LEGAL NOTES

COPYRIGHTS

All Material Copyrighted(c) 2001-2003 by the General Electric Company, All rights reserved.

REVISION HISTORY

Revision	Date	Reason for change
Rev 0	11/10/2003	501WB.3_H5.2M3

Page Left Intentionally Blank

TABLE OF CONTENTS

Legal Notes.....	5
Table of Contents	7
Features LightSpeed RT.....	11
SOFTWARE VERSION.....	11
GANTRY.....	11
SCAN CONTROL ARCHITECTURE CHANGES	11
OPERATOR CONSOLE/GLOBAL RECON ENGINE	12
HIPAA.....	12
CHANGES TO SHUTDOWN	12
SMARTPREP	12
REPEAT SERIES	13
LOW SIGNAL CORRECTION	13
FILM COMPOSER DEPOSIT TIME	13
COLOR IMAGE SUPPORT	13
AUTO TRANSFER.....	13
ANONYMOUS SCAN DATA SAVE	14
WIDEVIEW RECON ALGORITHM	14
EXTENDED CT NUMBER RANGE.....	15
HELICAL TILT	15
THIN TWIN HELICAL.....	15
HELICAL SCAN TYPE	16
CINE SCAN TYPE.....	16
AUTOMATIC EXPOSURE CONTROL – AUTOMA.....	16
New Precautions - 501WB.3_H5.2M3	19
SYSTEM.....	19
HIPAA LOG IN	19
SCAN.....	19
TUBE WARM UP/FAST CAL	21
PATIENT SCHEDULE	21
PERFORMED PROCEDURE STEP (PART OF CONNECT PRO OPTION)	21
PROTOCOL MANAGEMENT	21
SMARTPREP (PURCHASED OPTION).....	22
RECONSTRUCTION	22
3000 IMAGE SERIES.....	22
RECON MANAGEMENT	23

DISPLAY (APPLIES TO BOTH IMAGE WORKS AND EXAM RX DISPLAY)	23
IMAGE WORKS DISPLAY	23
EXAM RX DISPLAY	24
VARIVIEWER	24
DIRECT3D - OPTION	25
VOLUME VIEWER (OPTION)	25
CTCOLON (OPTION)	25
ADVANCED VESSEL ANALYSIS (OPTION)	25
REFORMAT	26
CT PERFUSION 2 ON OC – PURCHASED OPTION	26
ADD SUBTRACT	26
ARCHIVE	26
FILMING	27
NETWORK	27
IMAGE MANAGEMENT	27
LEARNING AND REFERENCE GUIDE	28
Legacy Precautions	31
SYSTEM	31
TUBE WARM UP/FAST CAL	32
SCAN	32
SMARTPREP (PURCHASED OPTION)	35
PATIENT SCHEDULE	35
CONNECTPRO (PURCHASED OPTION)	35
PERFORMED PROCEDURE STEP (PART OF CONNECT PRO OPTION)	35
RECONSTRUCTION	36
RETRO RECON	36
RECON MANAGEMENT	37
CT PERFUSION2 ON OC – PURCHASED OPTION	37
AUTO TRANSFER	37
DISPLAY	37
EXAM RX DISPLAY	38
IMAGE WORKS DISPLAY	39
EDIT PATIENT DATA	40
REFORMAT	40
FILMING	40
ARCHIVE	42
NETWORK	43

IMAGE MANAGEMENT	44
------------------------	----

FEATURES LIGHTSPEED RT

SOFTWARE VERSION

501WB.3_H5.2M3

This is revision 0 release of software for the LightSpeed RT with 80cm bore.

GANTRY

Gantry – has been redesigned to support an 80 cm opening with expanded DFOV to 65cm using WideView recon mode. The gantry has new tube/detector geometry. Focal spot to Isocenter is 60.6cm; focal spot to detector is 105.87cm.

The alignment light positions have been redesigned for improved positioning for oncology patient imaging.

Detector/DAS configuration is a 4-slice MDAS type supporting acquisitions in axial, cine and helical modes with 0.75:1 and 1.5:1 helical pitches.

Rotation times of 1, 2, 3, and 4 sec for Axial and 1 sec for Cine and Helical modes.

SCAN CONTROL ARCHITECTURE CHANGES

Table and Cradle - are always referenced. No need to reference the table after a system reboot.

After a shutdown of power to the Gantry, on start-up the E-Stop reset **MUST** be pressed for the gantry control panel to reset the hardware.

The following sequence should be followed:

1. Turn on A1 power
2. Press reset button on gantry controls
3. Turn on OC console and the system will complete the start up.

During start-up of the system, the sequence of the hardware is as follows:

- o Gantry display is active even after system software has been shutdown,
- o Gantry display resets after applications begin to start up,
 - Lights will sequence through the display,
 - Scan hardware is reset when interference display stops flashing,
 - If alignment light is still flashing after applications have completely started then hardware is not up and a scan hardware reset will need to be done.

Activation on one of the collision tape switches on table will stop down table motion. A reset is no longer required. Move the table away from the obstruction, clear the obstruction and then the table may be lowered.

E-Stop removes power from the table.

- o Select reset, will be ~15-20 second delay for the reset to complete before table is operational.

OPERATOR CONSOLE/GLOBAL RECON ENGINE

Xtream™ Technology - The operator console is now uses a PC based computer system running a Linux based operating system. The system supports a gigabyte network port. The system supports 3000image series and has on-system storage for 250,000 images.

GRE Recon – The Global Reconstruction Engine hardware provides faster reconstruction times of up to 6 images per sec. Changes in the recon and image have support a time to first image of 2-3 seconds. The faster times enhance the capabilities of SmartPrep to provide more precise timing in the delivery of contrast material for a wide variety of examinations.

Applications on OC – Common user interface with Advantage Windows for improved productivity. The OC now supports an expanded list of Advanced Application packages such as Volume Viewer, Advanced Vessel Analysis, CT Perfusion2, DentaScan.


HIPAA

HIPAA (Health Insurance Portability and Accountability Act) -The system provides the option to require users to log in and enter a password to use the system. This capability has been added to enable your institution to comply with HIPAA. Users may be assigned to groups of users and permissions may be set for each user. Permissions may also be set to restrict access to the Protocol Management feature. Your field engineer configures log in capability. Login is selected from the shut down button.

CHANGES TO SHUTDOWN

Shutdown/Restart/Login - When shutdown is selected, you will be presented with 2 choices Restart or Power down. Restart shuts down, and then automatically restarts the system. Power down brings the system down to the power off prompt. When shutdown is selected, you must cycle power to the operator console to reboot the system.

If the HIPAA login feature is enabled, you will be presented with 3 choices Logout, Restart or Shutdown.

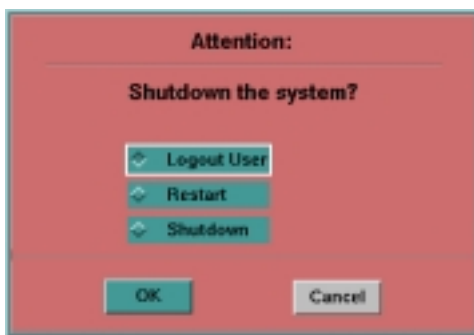


SMARTPREP

SmartPrep Monitor ISD - With faster image-to-image times the minimum Monitor ISD has been decreased to 1 sec. This provides faster monitoring of contrast enhancement. The number of Monitor scans has been increased to 40.

REPEAT SERIES

Repeat Series – has additional functionality added in this release of software. Repeat Series now allows you to choose any series that has been scanned. When more than 1 series has been scanned, a list of all scanned series will be displayed. Click on the series that you wish to repeat. If only one series has been scanned, Repeat Series will not display the Repeat Series pop-up. The state of Show Localizer on or off is remembered by the system. If Show localizer state is on, it will be automatically displayed when you return to the ViewEdit screen.

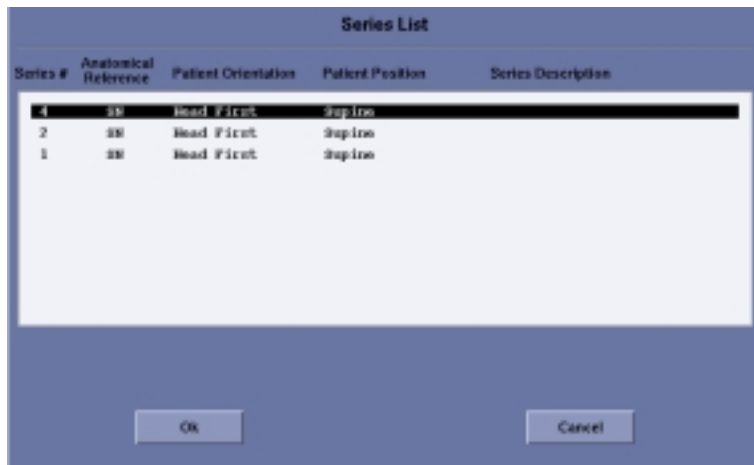

LOW SIGNAL CORRECTION

Advanced Artifact Reduction (AAR) - reconstruction correction algorithm has been updated for this release. A reduction in streak artifacts due to low signal condition through shoulders and hips should be seen.

FILM COMPOSER DEPOSIT TIME

Auto Film Composer - The time to deposit an image on the Auto film composer has been reduced to 350ms. This will allow Auto Film to improve productivity in conjunction with the faster reconstruction times.

COLOR IMAGE SUPPORT

Color Image support - is provided for the sending and receiving of color images created in the Advanced Applications.


AUTO TRANSFER

Auto Transfer - by image has been updated to transfer in larger groups of images. This provides more efficient use of system resources while Auto Transferring by image. Auto transfer by image will transfer the first 50 images in 2 groups of 25, after that images will be transferred in groups of 50.

ANONYMOUS SCAN DATA SAVE

Scan Data Save - The system now provides the ability to save scan data with patient identifying information removed. This has been added to address HIPAA concerns. Raw data saved anonymously and then restored will reconstruct with patient identifying information removed. The Patient Name and ID will be anonymized.


- Select - Recon Management
- Select - Save Scan Data
- Select - Save Selected Annon. Data

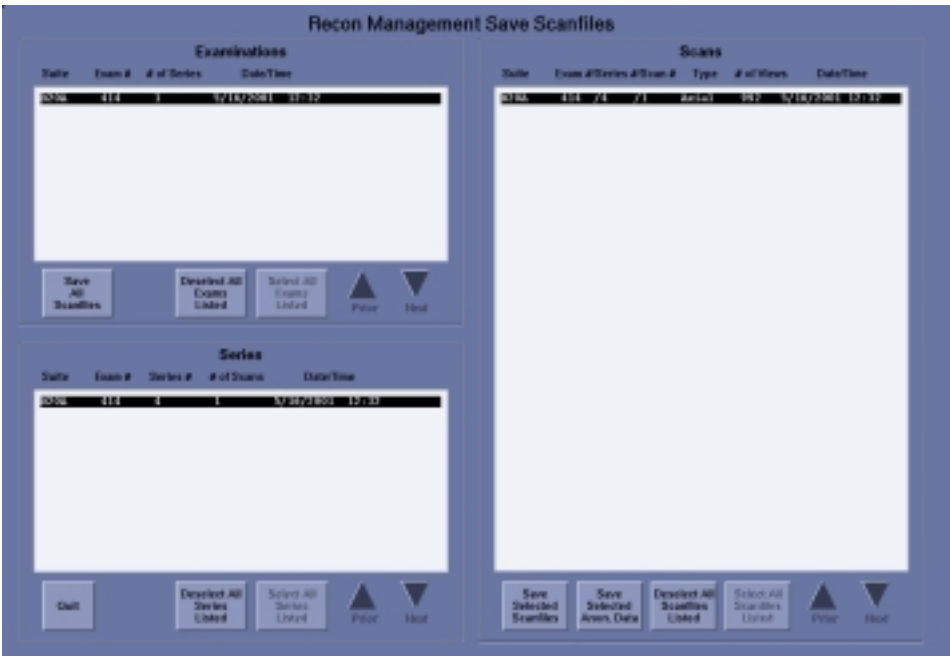

WIDEVIEW RECON ALGORITHM

WideView Recon Algorithm— supports reconstruction of images with a 50.1 to 65cm DFOV when acquired with the Large SFOV.

The Edge recon algorithm is not compatible with WideView Recon.

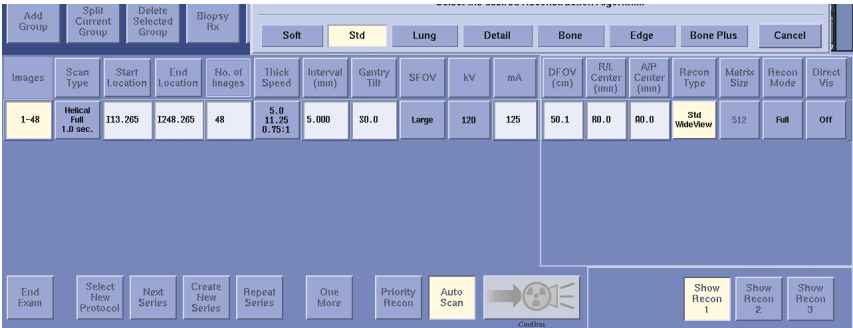

EXTENDED CT NUMBER RANGE

CT Number Range – expanded to support +/- 31743 pixel values in the image. The expanded CT Number range will help in review of images where metal implant such as surgical clips, dental work or joint prosthesis exist.

HELICAL TILT

Helical Tilt – provides tilt capability for Helical scan type. The system automatically adjusts the table speed to maintain prescribed interval based on the degree of tilt prescribed. Images are annotated with the prescribed slice thickness and scan location. The table speed shown is the adjusted table speed depending on degree of tilt. Helical Tilt is not allowed for slice thicknesses of 3.75mm or greater with PED Head and Head Scan-Field-of-View (SFOV).

Images	Scan Type	Start Location	End Location	No. of Images	Thick Speed	Interval (mm)	Gantry Tilt	SFOV	kV	mA	Total Exposure Time
1-33	Helical Full 1.0 sec.	126.540	353.460	33	2.5 3.75 0.75:1	2.500	39.0	Head	120	200	23.8

THIN TWIN HELICAL

Thin Twin Helical – acquisition mode supports the implementation of sub-millimeter (0.625) imaging for the helical scan type. All recon types are supported with Thin Twin helical mode. The pitch is 1:1 with fixed table speed of 1.25mm per rotation. Thin Twin acquisitions are limited to the small focal spot. Gantry tilt is allowed. Thin Twin acquisitions cannot be combined into a thicker slice thickness in Retro Recon. ThinTwin is not available in Biopsy mode. Thin Twin is compatible with all SFOVs and Direct3D.

Select the desired Image Thickness

Detector Rows		Detector Configuration:	
2	4	2 x 0.625	
Helical Thickness (mm)		Beam Collimation:	
0.625	1.25	1.25mm	
2.5	3.75		
5.0	7.5		
10.0			
Pitch		Coverage Time:	
1:1		65.3 sec.	
Speed (mm/rot)			
1.25			
OK		Cancel	

HELICAL SCAN TYPE

Helical Segments – are supported retrospectively in Retro Recon in this release of software.

CINE SCAN TYPE

Cine Segments – are supported prospectively in Scan Type and retrospectively in Retro Recon in this release of software.


AUTOMATIC EXPOSURE CONTROL – AUTOMA

AutomA – automatically adjusts the mA used per rotation along the z axis of the patient. The amount of adjustment is determined from a single scout image, the last scout scanned, and is calculated based on the Noise Index entered. Each Noise Index is optimized to the slice thickness and the anatomical region the protocol was selected from. Auto mA can help to reduce dose to the patient by optimizing the mA used during a scan sequence.

There are three tables whose Noise Index values are 15% different from one another. Table 2 is loaded at software installation.

The mA button is now a pop-up screen, which support both AutomA and Manual mA.

mA Control pop-up

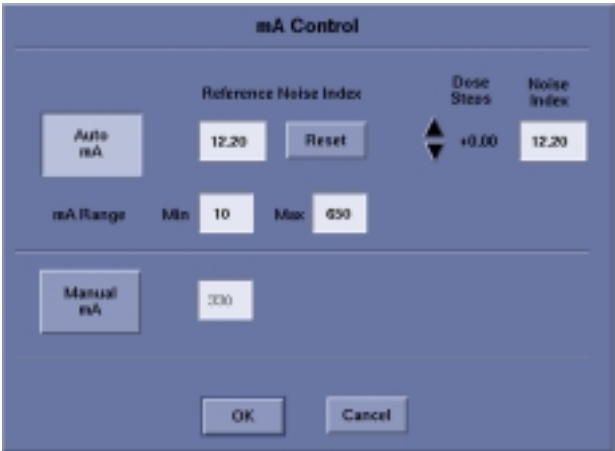

Manual mA	The Manual mA button is used to set up the system so mA values can be entered manually in the text field for mA.
Auto mA	The AutomA button turns AutomA on. In Protocol Management, the Reference Noise Index, Min/Max mA Range needs to be defined for the protocol. In ViewEdit during scan prescription, the Noise Index or the Dose Steps can be used to modify the Noise Index for the current acquisition.
Reference Noise Index	Default or baseline Noise Index for the given protocol. Any changes to Dose Steps, Slice Thickness or Noise index are referenced to this value. This value can only be prescribed while in protocol management.

Noise Index	The noise level required for the study. As the Noise index increases the required mA decreases and image noise increases.
Dose Steps	Adjusts Noise Index by steps of 5%. Dose steps can be increased or decreased + values decrease image noise thus increasing required mA. Minus values increase Noise index, thus decreasing required mA. A Dose step value of 0 indicates that the prescribed Noise index is equal to the Reference Noise index for the protocol.
Reset	Resets the Reference Noise Index to the GE Target Noise Index Default for the anatomical area and slice thickness chosen in the protocol

Noise Index will be annotated on the image as NI:(value).

The average mA for the all views which are used to create the image will annotated as the mA value.

You can find more information on AutomA in the Building Protocols chapter in the Learning and Reference Guide provided with your system

This page left intentionally blank

NEW PRECAUTIONS - 501WB.3_H5.2M3

SYSTEM

- The scan hardware may fail to reset on reboot. Select Service Desktop, System resets, Scan, Run to reset the scan hardware. If the reset still fails, then shutdown and reboot the system to recover.
- The gantry display may show gantry tilt at S0.5 or I0.5 even though the gantry was at zero when the gantry was powered down. You may have to return the gantry to 0.
- The system may not show typed information on the New Patient or ViewEdit screen. The system will update after a short time.
- The system may be slow to respond to keyboard entry and mouse clicks. If the system fails to respond stop and let the events catch up.
- The dynaplan screen may fail to update if Priority Recon is selected during scanning
- Pressing the Space Bar after selecting all fields in New Patient locks the patient information fields and all the patient information fields become insensitive. Only the delete and backspace key will be active. Press the delete key to recover from this situation.
- Entering a \ and pressing the Enter key in the New Patient can cause the fields to become insensitive.
- Under heavy system operation, the message "Unable to install image in the database" may be seen. Recon for the series will take longer due to unsuccessful install of images to the database.
- A pink attention box may not go away even after OK is selected. If after 15 min, the attention box has not gone away, go to the Image Works desktop and select Unix Shell from the Tool Bar in the upper right corner of the monitor. Turn Cap Lock off. Place the mouse in the Unix shell window, and enter the following command: "slay att". The system will echo a response for the attention box that remains open, enter a Y for yes. This will close the attention box.

An open attention box will cause the system to become very slow until the message box is slayed.

- If "err" is seen on the gantry control panel after a reboot or system reset, perform a system reset.
- Coronal laser light accuracy is +/- 3mm at 1 degree. Improvements are being investigated for future release.

HIPAA LOG IN

- To move between the user name and password field press the tab key twice. Selecting enter may reset the user password inadvertently resulting in an invalid password.

SCAN

- Show localizer does not show the Scan Field of View when prescribing scans. When scanning head studies, verify that you have your patient centered at ISO center to assure that the entire head is included in the SFOV.
- Show Localizer may fail to display. If toggling the Show Localizer button does not display the scout, then display the scout in a free viewport and use Report Cursor to find your scan locations and A-P and R-L centers, manually enter these on the ViewEdit screen.

Scan cont.

- Show Localizer may remove the localizer lines when page down or page up is selected and will not advance to the next scout. To recover from this, toggle Show Localizer off / on. If Scout Image 2 is displayed use Page Up, to go to scout image 1 use Page Down.
- The Show Localizer button may get in a stuck state and not display a scout image if the last scout scanned has an azimuth other than 0, 90, 180, 270. Manually display a scout in a free viewport and use report cursor to find your scan locations and then manually type in the Start and End location, RL and AP Centers as needed.
- For DFOV of greater than 49cm, the graphic Rx handles for tilt and DFOV width may not be able to seen or selected. Adjust the DFOV, A/P and R/L centers manually in the Recon Tab on the ViewEdit screen.
- The pause button on the screen in Tube warm-up does not function. Use the stop scan button on the SCIM.
- The table will take 10mm to stop if the table is running at 100mm/sec when a patient weighing 450 lbs. is on the table.
- You may be prompted to run DAS gain cal prior to starting Fast Cal.
- The system may not allow changing of the landmark even though no scans have been taken in the series. A landmark that is set at the edge of the scannable range most often causes this. To reset the landmark, choose Select New Protocol.
- A scan abort may occur, resume to continue. This can be caused by the following:
 - Taxi link violation due to forward error correction
 - DIP waiting to receive buffer of views
 - Generator Errors
- The system does not verify patient weight against the selected Pediatric protocol when the protocol number is typed-in. Always verify that you are entering the correct protocol for the weight category of the patient.
- If Direct3D or VariViewer are enabled in DirectVIS under the Recon tab, Add Group will display some fields as insensitive as it will be combined with the current Direct3D or VariViewer session. Turn Direct3D or VariViewer off for the added group if you no longer want it to be included as part of the Direct3D or VariViewer session. This will allow changes to any of the acquisition parameters.
- Occasionally an image may be reconstructed such that it is displayed in a rotated condition. Go to Retro Recon and reconstruct the image again.
- Scan files restored and reconstructed from prior release of software may reconstruct the images in a rotated position. This is because the start angle was not being recorded and now it is recorded as part of the scan file information.
- Selection of the Edge reconstruction algorithm with WideView recon mode can cause recon to fail. Do not select Edge for DFOV greater than 50cm.
- Truncation artifacts may be seen at the edge of image data that violates the 50cm SFOV. These truncation artifacts do not go away even when WideView recon is enabled.
- If large off-center field of view center is targeted, the CT numbers on the side that is outside the DFOV may be seen a -2000. Do not use small DFOV with large off-center targets.
- If there is a large gap between tissues at the edge of the SFOV for a WideView reconstruction, the IQ may be degraded in this area. Place a saline bag in the area to eliminate the air gap.
- Pressing the Prescribe Tilt button on the OC again after the gantry has reach the specified tilt angle will cause the gantry to move back towards zero. Once the gantry is at the correct tilt angle do not continue to press the Prescribed Tilt button.

Scan cont.

- If the table is set to the lowest position of 264 from isocenter and helical scan is prescribed, the next group may not be able to be scanned because the table will no longer be in a scannable range. Make sure the table is at 263.5 or above to prevent this situation.
- If a large tilt angle is required for a helical scan, it may be better to use axial scan type, or if helical is desired, scan with no tilt and reconstruct in the tilted plane in reformat for optimum image quality.
- GE protocol 22.9 has filming turned on. Go to Filming tab and turn AutoFilm off if you do not wish to film. GE protocol 25.15 in Chest is corrupt and cannot be copied. Copy protocol 26.15 in Abdomen to Chest and adjust mA accordingly.
- If the last scan before selecting Biospy mode is a Cine scan type, system will allow the Cine scan type. Only axial and helical scan types should be used for biopsy mode.
- The Axial dose reported in the Technical Reference Manual is incorrect. The body dose information is reported as 5.5 mGy/100 mAs center and should be 4.2 and is reported as 10.0mGy/100mAs center and should be 8.3.

TUBE WARM UP/FAST CAL

- At the beginning of Fast Cal the system will do a balance check. This will take approximately 2-3 min and then Fast Cal will begin.
- Fast Cal may give the error that Collimator Cal failed. Resume Fast Cal and contact GE Service.

PATIENT SCHEDULE

- The Patient Schedule button may not display the Work List from the HIS/RIS server if the network is slow. Try again.

PERFORMED PROCEDURE STEP (PART OF CONNECT PRO OPTION)

- The PPS server may stop and fail to send the request. You may have to complete multiple times to get the series to update to the complete state.
- Add Sub does not update images for PPS correctly
- Perfusion2 will not launch if images contain PPS information.

PROTOCOL MANAGEMENT

- The Dose information in protocols used from Most Recent does not reflect the dose actually used for the protocol. If the protocol contains manual mA, retype the mA value to update the dose display. If Auto mA is used, the dose will not be reflected in Protocol Management, but will be updated at scan time for the patient being scanned based off a valid scout image.
- There is no default protocol for Pediatric areas Neck, Upper Extremity, Chest, Abdomen, Spine, Pelvis, and Lower Extremity.
- If the weight of a child is at the cross over point of a weight category, due to rounding errors the correct weight category might not be selected. Please check the weight based category selected against the label and the patient's weight.

SMARTPREP (PURCHASED OPTION)

- The clock countdown in the lower right viewport of the SmartPrep display starts counting at the press of Monitor Phase. The countdown should begin at the press of Start Scan. To avoid a large offset in the time from start of injection, do not select Monitor Phase until the injector is armed and you are ready to begin.
- The transition from Monitor phase to Scan Phase will be longer than the delay prescribed in SmartPrep set up. Actual minimum delay is approximately 6-8 seconds.
- SmartPrep may fail to confirm to take the baseline image. Try confirming again, if this does not clear the problem then a system reboot will be needed to get SmartPrep working again.
- The SmartPrep baseline image may fail to display even though it is recon'd and in the database. Go back to ViewEdit and rescan the baseline image to continue. The Monitor phase image and graph information will not update if the table height has changed from the Baseline to Monitor phase. Do not change the table height after the baseline scan has been taken to avoid this.
- When transitioning to the Scan Phase from Monitor phase if the XR Auto View layout is selected the upper left viewport will be black. Move your mouse cursor over the viewport to redisplay the image.

RECONSTRUCTION

- Streak artifacts may be seen in images when a tube spit has occurred.
- Image recon may fail and then shutdown. Select Recon Management, select Pause queue, Restart queue, if this does not restart recon then a Restart of the system will be needed.
- Cardiac recon modes of Burst 2 or Burst 4 cannot be selected via column edit. Please select the individual row to edit recon mode selection.
- The dynaplan screen may fail to update if Priority Recon is selected during scanning
- Retro recon may not be able to get the same image locations as prospective recon. This is due to rounding in the start and end location. To avoid this mismatch prescribe start and end locations that are even numbers.
- Retro recon list service exams from Fast Cal. Do not attempt to recon this data, the SRU will shut down. These exams are listed with exam numbers that begin with 50000.
- Anatomy scanned off center from ISO center such as extremities, may have loss of resolution when data is acquired helically. Always make sure to center to anatomy of interest as close to ISO center as possible.
- Cine acquisitions may have retries or image suspension. Unsuspend images in recon management.

3000 IMAGE SERIES

Display

- Paging and image review with next prior may slow when reviewing large series.

Archive

- An image series of 3000 images can only fit on one side of a 2.3 gb MOD.
- Auto Store will not save the exam if the media size attached is 1.2 gb.

3000 Image Series cont.

- If a 1.2 gb media is used for image archive and the image series is greater than 1500 images archive will need to be queued by image. Select a range of images 1500 or less when archiving to 1.2 gb media.
- If greater than 2000 images are queued to save by image you may see the browser disappear and reappear. While this is occurring it will be impossible to access the browser. To avoid, this queue less images at a time when saving by image.
- Images may fail to save if saving large groups of images. Verify that all desired images are saved before deleting them from the system disk.

Network

- When queuing a large number of entries, there may be delays when accessing the network queue. Typically this will be seen when 200 or more entries exist.

Dentascan/Add Subtract

- These applications can only load 1000 images. Highlight the desired range of 1000 images by selecting the first image desired, hold the shift key down and highlight the last image desired then select the desired application

RECON MANAGEMENT

- Scan File save is not complete even though the Saved File message has been posted. The problem is due to the time it is taking to unmount the DVD. You may have to wait 5-7 minutes for the DVD to unmount. The write light on the DVD will be green when the DVD is unmounted
- Scan Files restored from a different system will be reconstructed using the product name of the system it is restored on. For example, restore scan files from a H1 SDAS QX/i or a H2 SDAS Plus to a LightSpeed¹⁶ and the product name on the images will be indicated as LightSpeed¹⁶ instead of LightSpeed QX/I or LightSpeed Plus. All other patient and parameter information is correct.

DISPLAY (APPLIES TO BOTH IMAGE WORKS AND EXAM RX DISPLAY)

- Plus recon mode is not annotated on the Series Text Page. Refer to the image annotation for the Plus recon mode annotation.
- An active user annotation graphic will not film the box or arrow that is displayed on the screen.
- An arrowhead cursor will be displayed instead of crosshair cursor in the AutoLink, Free viewports on Exam Display desktop and on the viewports in the Viewer on Image Works desktop.

IMAGE WORKS DISPLAY

- The primary focus may change if the space bar is selected multiple times when entering accelerator line commands. Verify that primary focus is on the image desired.
- Entering Quit in the accelerator line in the Viewer will close the Viewer; however the system will post the error message "Viewer failed to start."

EXAM RX DISPLAY

- The Trackball may fail to initialize properly at system reboot. Try rebooting the system again to recover
- The Icon for paging/WL and Auto Link is not visible
- Display windows may display on the left screen instead of the right. This will most often occur when switching desktops rapidly. A system reboot will be needed to recover.
- Pressing page up 2 times causes the ww and wl to change on the scout on it's own. Readjust the ww and wl to continue.
- Images reconstructed with Priority recon are not available in the Auto Link viewport. Redisplay the series to access these images.
- The first time the trackball is used for paging after a system reboot window level will be changed instead of paging the images even though the P icon is displayed. Click once to switch to W mode then click again to place the trackball in paging mode.
- The AutoView or Auto Film viewport may display a very magnified image. This image may be filmed in Auto Film. Review your films and verify that all filmed images look correct.
- Viewport format may not be able to be changed and requires switching layouts to correct. This occurs if F3 or Shift F3 is used for filming and next image is used. The viewport format is now locked. You will have to switch to another review or Auto View layout to correct.
- MIROI may fail to produce a graph if next prior is selected while doing next prior. Redisplay the series and begin again.

VARIVIEWER

- Batch film may fail after a system reboot if the format selected is different than the default film format on reboot. An error dialog will be posted informing you of the error. You may have to select filming multiple times to successfully send the request.
- If the scout is taken after the series used for VariViewer operation, the system may be unable to find the scout for graphic batch prescription. Manually enter the range for the batch prescription.
- Any zoom applied to a Variviewer image is not maintained as you page through the images.
- VariViewer session may not be able to be started if there are multiple Pause scans during the acquisition of data.
- VariViewer is not compatible with a gantry tilt.
- Switching back from a 768 viewport is VariViewer causes the VariViewer viewport to be forward instead of the Auto View viewport. Use the page turner to toggle back to the auto view Viewport.

DIRECT3D - OPTION

- Batch film may fail after a system reboot if the format selected is different than the default film format on reboot. An error dialog will be posted informing you of the error. You may have to select filming multiple times to successfully send the request.
- 768 image size may display split in 2 if resize is selected with the Auto View viewport forward. Always have the Direct3D viewport forward when selecting resize. To recover from this problem select resize twice.
- If you resize the viewport in Direct3D, the image may not initially be displayed. User the Page Turner in the lower left corner of the viewport to toggle to the time.
- The Direct3D model may not initially build properly. Enter Interactive Review and re-render the model and then the model will build correctly.
- SaveState for Direct3D does not load into Volume Viewer on the OC. It does work correctly when sent to the AW.

VOLUME VIEWER (OPTION)

- Volume Viewer may display on the ExamRx desktop. This will occur if you switch back to the ExamRx desktop quickly after launching Volume Viewer
- Volume Viewer may hang. This occurs due to an error dialog that is hidden. To avoid this, don't switch desktops until you start to see the model build in Volume Viewer. This may also occur when changing anatomical areas. The Volume Viewer software incorrectly may think the data set is not valid.
- A 3D model may be missing data when it is initially displayed. Slightly rotate the model to redisplay all of the 3D data.
- The system may get extremely slow when saving a 3D or VR model. Please be patient for the system to complete saving the model.

CTCOLON (OPTION)

- Rotate and translate is missing in the CTC and Virtual Dissection protocols
- The system may hang and display 3 film composers with no borders, the composers cannot be moved and no selections can be made except for desktop switching. A system reboot will be needed to recover.

ADVANCED VESSEL ANALYSIS (OPTION)

- AVA defaults to inferior view instead of Anterior after processing AVA. Use the anatomical buttons to change the view to Anterior
- AVA lumen view may show wavy lines while rotating. When the rotation is stopped the display will update.
- The system may hang and display 3 film composers with no borders, the composers cannot be moved and no selections can be made except for desktop switching. A system reboot will be needed to recover.

REFORMAT

- Reformat will only allow filming in Batch to the first 30 cameras in the list. If another camera is desired save the batch images to disk and use Print Series to film the images.
- Images saved from batch mode may appear to have jagged edges or not have optimal Image Quality. Manually save the images to improve image quality.
- Large image series may fail to load when using Reformat detail. This can be caused by small spacing interval between images or image sets greater than 1000 images. Select a sub set of images if this occurs or use Reformat Standard to build the model.
- A series that has more than 1 group with 0.625mm slice thickness will fail to load for reformat or will fail in Volume rendering on the AW. This is due to a mismatch in the spacing between the groups. To avoid this have only one scan group if possible, otherwise selectively highlight the images from the first group for Volume Analysis.

CT PERFUSION 2 ON OC – PURCHASED OPTION

- Computation of perfusion maps and automatic detection of the artery or vein may not function when the data displayed in Perfusion is initially black and the WW and WL values need to be set to large negative values (WW 40, WL -1020) to make the images viewable. The threshold should be set to the maximum negative value for air to assure that automatic detection of the artery and vein and that Perfusion maps are displayed correctly.
- Perfusion will not launch if the images contain PPS information
- Saving TIFF files of the graph and images is not supported because there is no floppy drive hardware available on the OC.
- Set the color ramp to gray scale for screen saves images if you are going to network the images to another station such as PACS to improve quality of the screen save.
- Initially as CT Perfusio2 is loading, you may see garbage flash along the left side of the screen.
- Processed data does not display if New Protocol is selected within Perfusion2. To avoid this always close Perfusion2 before selecting a new protocol.

ADD SUBTRACT

- Comb images will not contain the e/s/i used for the added image.

ARCHIVE

- Even though the Sony MOD drive states that it is a 5.3 gb drive. The archive software only supports a maximum size of 2.3 gb. 2.3 gb 512 bytes per sector is the preferred archive media size to utilize on the system.
- Archive may fail and slow the system when more than 200 images are queued by image. Try to archive by series if possible to avoid this.

Archive continued

- The system may report that the media is full even though the media has just been labeled when Save by image is selected. Try selecting a smaller range of images to avoid this.

FILMING

- Image setting such as WW WL and flip rotate, zoom, roam are not maintained across all images when imitating F4 print series from a MID Viewport format. Use a 512 size viewport to initiate F4 Print Series.
- Auto Film may stop and the Continue Same Sheet button will be displayed. The Continue Same Sheet button will be active, but when selected will not restart filming. Use Manual Film to complete filming the exam and reboot the system before starting the next exam
- If the F1 key is used to film the Text Page ROI, Exam Text Page and Series Text Page, the text pages will be filmed to the film composer last selected in the pop-up for text page. Use the filming selections in the text page screens to film the information displayed.
- The message "Film formatting in progress. Please retry." May be seen when filming to the Manual Film composer when the system is busy with other simultaneous operations such as recon, network, auto film and scanning. Retry manual filming when the system is less busy. A system reboot may be required to recover.
- Image annotation may be small or blurry on filmed images when the first image is a scout that is over 600mm. Place this image on the last space of a film to avoid this.

NETWORK

- The system may create core files that may fill up disk space on the system. Networking color save images from Volume Viewer or CardiacIQ causes these core files. To avoid this save the color images to MOD, remove the color images from the system then restore the color images and then network the images.
- Color Saved images from Volume Viewer or CardiacIQ may fail to transfer. To avoid this save the color images to MOD, remove the color images from the system then restore the color images and then network the images.
- Images that have been edited using Edit Patient networked to a RadStore system will fail to archive to the RadStore DVD. The images must be archived on the CT system MOD.
- The system no longer supports Advantage Net protocol. When sending images DICOM send should be utilized.

IMAGE MANAGEMENT

- Removal of images should be done with the system in an idle state. Failing to remove images at idle times can cause image install timeouts that can then cause recon to shutdown and can cause the browser to fail to respond. The user may also see a blank list in the browser. The system should not be shutdown until the image space has been updated fully from the remove. The length of time to synchronize the database is proportional to the number of images being removed. This will assure that all images are reconstructed, displayed and installed into the database without error and that image space is update correctly.

LEARNING AND REFERENCE GUIDE

- After a desktop switch back to Learning Solutions desktop the guide returns to the home page instead of returning to the page they were on when the switch occurred. The document is slightly behind the home page. Left click and drag on the home page border to move the home page so you may access the Learning and Reference Guide pdf file.

This page left intentionally blank

This page left intentionally blank

LEGACY PRECAUTIONS

SYSTEM

- Run Shutdown and Restart daily.
- The system may fail to start if system disk space needs to be recovered from unused system files. The following message will be posted in the Disk Management shell window:

```
"The purpose of storelog is to recover system disk space by archiving and then removing core, log, and data files that have been saved for their system troubleshooting diagnostic value. Now they may be taking up too much space for the system to run properly. Removing these "system log" files does not add image space, but should allow the application to start."
```
- You are then directed to answer the following questions (be patient for the messages to display):

Place the mouse cursor inside the black Disk Management window

Q: Do you want to save system log files to removable media?

Type n and hit the enter key

Q: Do you want to remove the system log files?

Type Y and hit the enter key.

You will then see a message "Please be patient, the system is now removing unwanted files." The system will then start normally.

- If the system fails to startup completely, select [Unix Shell] from the tool chest menu in the upper right of the Scan Monitor place your cursor in the blue shell and type st.
- If mouse selections fail, select the [Escape] key to clear backlogged requests.
- If you can no longer type in a field, move the cursor to the field and try clicking the middle mouse to restore functionality. If this does not work, then a system shutdown will be required.
- The cursor may not move between the left and right screen. Wait a few seconds without moving the cursor and try again.
- To install a SMPTE pattern, BRH or Quality assurance images for viewing, select the following. Once installed the images will be available for viewing from the Image Works browser or List Select in Exam RX. Both will be listed as Exam 1000, the patient name will reflect if it is a QA image or a SMPTE pattern.
 - Select Service Desktop
 - Select [Diagnostics]
 - Select [Display Processing]
 - Select [Install SMPTE Image]
- In general, wait for a screen transition to take place before making another selection.
- If the console becomes unresponsive for 2 minutes or more, shutdown the system using the pink shutdown button and restart the system. If you cannot select the pink shutdown button, turn off the console power switch, wait 10 seconds, then turn the console power switch back on. The system should come up normally.

System cont.

- The system will lock up if the mouse becomes hung on the left monitor. The only known recovery is to cycle power to the console.
- You may see a pink warning dialog posted indicating the system has low disk space. This is due to a partition on the system disk getting too full. This can be caused by archive, network or film queues containing too many entries. Make sure that entries are processed. If the queue is paused, resume it or delete entries as needed.
- Protocol Management may not be selectable. Check that Retro Recon is selectable. If it is not; then the system will need to be rebooted.
- Pop-up screens and menus may appear on the wrong monitor or may be displayed split between the two monitors.
- Windows such as the Film Composer that normally are restrained to the right or left monitor can be dragged to the other monitor. Be careful not to get a Film Composer hidden behind some screen if you move it to the other monitor.
- If a blank viewport is seen where image annotation is displayed and W/L is interactive but no image is seen, shutdown the system, recycle power and reboot the system to clear this display issue.
- If corrupted images are seen after applying an image filter such as enhance, smooth, lung or gray scale, try rebooting the system and see if the image can be displayed correctly. If not contact your GE service representative.

TUBE WARM UP/FAST CAL

- Run Tube Warm Up after 2 hours of non-use. It is best to warm the tube as close as possible to the time the next patient will be scanned. Don't warm the x-ray tube and then wait 30 minutes without scanning.
- Always check that the beam is clear when doing Fast Cal. In Fast Cal, the system checks cleanliness of the Mylar Window. If the system suspects the Mylar Window may be dirty to the point that it may cause a beam obstruction, a pop-up message will be displayed. The operator will be instructed to clean the Mylar Window and select [Retry].
- Fast Cal must be run once every 24 hours.
- Complete all portions of Fast Cal, Warm Up 1, Warm Up 2, Gen Cal, Clever Gain and Fast Cal. This assures that the Air calibration and generator calibrations are up to date on the system. Scan aborts may occur during Axial or Helical scanning. Always be aware of the scan progress during an Exam and select Resume as soon as it is posted to continue.

SCAN

- Scan aborts may occur during Axial or Helical scanning. Always be aware of the scan progress during an Exam and select Resume as soon as it is posted to continue.
- Scan may fail to confirm posting a message that not enough image space exists, even though the image space shown in the Feature Status Area indicates there is enough space. This is due to the fact that images are stored on the system disk in more than one partition. Remove consecutive exams to free up image space for confirm to proceed.

Scan cont.

- To terminate an Insite connection on your system. Select [New Patient] a message will be posted informing you that the scan hardware resource is not available. Wait 3 minutes and select [New Patient] again to begin scanning.
 - In general, if a scan fails and a [Resume] is posted, press [Resume] to continue. Try [Resume] again if the first [Resume] fails. If a failure still occurs, reset the scanning hardware through System Resets in the Service desktop. If scan still fails to restart, shutdown and Restart the system.
 - Auto Voice may fail to function, especially during system simultaneity. Make sure that you can hear the Auto Voice to recognize if Auto Voice has quit, manually breathe the patient when this occurs.
 - The Show Localizer scout image may fail to display or images may fail to recon and may become suspended if scanning is started while remove images is in progress. To avoid this **don't** remove images while an exam is in progress or scanning is active. Remove images when the system is idle.
 - The Dynaplan screen may report incorrect status (i.e. a scan is removed before it is actually scanned) when stop scan is selected. The screen will reset correctly after Resume or back to ViewEdit is selected.
 - Scan groups will not be contiguous if you switch from 2 or more Helical scan groups column or row edit and the scan type is changed to Axial. If this change is made, then verify the Start and End locations of each of your scan groups and adjust if needed.
 - The Cine time between images may change when other parameters are modified within a protocol. Check the Cine time between images prior to confirming the scan to verify it is the value you wish.
 - Tilt handles will not be visible on the screen if the DFOV is larger than 48. Use a DFOV smaller than 48 to assure that the tilt handles are visible.
 - Show Localizer may fail to display the scout image if the Next Series is selected before the scout images are reconstructed. Wait until the scout images are reconstructed before selecting Next Series. If Next Series has already been selected, toggle Show Localizer off and then on again.
 - Quickly changing desktops when Show Localizer has been selected for the next series can cause a partial display of the scout in the Graphic Rx window to occur.
 - If Page Up or Page Down fails to change the display Show Localizer:
 - Try the following in order:
 - 1. Turn Show Localizer OFF and then ON.
 - 2. 2. Move the cursor off the image window and back on to the image to refresh the screen.
- If both of these steps fail, then verify you have a valid scout for Show Localizer.
- If an R, L or A, P value is more than half of the Display Field of View in mm. The image annotation will not annotate R and L or A and P on the images. For example if the DFOV is 10cm and the R value is 56mm then the image annotation will show R R. Use a R-L or A-P value which is less than half of the DFOV in mm to avoid this.
 - A protocol may fail to display. The following message will be posted to the message line on the lower center of the scan monitor "Can't read selected protocol, please choose another protocol". If this occurs the protocol is corrupt and will need to be rebuilt. Delete the protocol and rebuild it.
 - The patient history field may be missing the last 9 characters when all 60 characters have been entered and the patient record is either selected as a completed exam or the Patient ID is entered and the system matches a completed exam.

Scan cont.

- Even weights entered in Patient Schedule may be rounded up to the next odd pound. This occurs because weights are stored as kilogram units then converted back to pounds.
- During the acquisition of large data sets (1000-1500), there may be Auto film and reconstruction performance issues. If either stops, check the queues and restart to continue.
- When selecting the Auto Voice field, the selected Auto Voice will be deselected. You must reselect the desired Auto Voice message before exiting the Auto Voice screen.
- When entering birth dates on the New Patient or Schedule Patient screens the following rules apply:
 - 2 digit years can only be entered if the birth year is in the current century for the 21st century - for the year 2000 forward.
 - Birth dates from any other century must be entered as 4 digits for example for the year 1899 –1999 all four digits must be entered.
 - Birth dates can only be entered for 150 years (Current year minus 150 years).
- If InSite is running the “Remote Safety” test, New Patient will not open. New Patient can be made available by either:
 1. Calling InSite and request they abort the test,
 - or,
 2. Go to the Service Desktop and select [Clean Up] option to can the test in progress.
- PMR images may fail to recon if disk space is low. To recon the images, remove images and restart the recon queue in Recon Management.
- You may get the message “Duplicate Scan Key”, if a Scout Scan is aborted within the first 20 to 75 mm of table travel or if a the prescribed length of the Scout is less than 75mm. In order to proceed, an End Exam must be done and new exam started.
- It may be possible to prescribe a recon range that is outside the recon 1 boundary. This may occur if the start location is dragged below the recon end location in recon 2 or 3.
- If a large number of scan files are reserved and there is a large recon queue, you may see the message: “Duplicate Scan Key – unable to allocate scan file.” There will be a resume at Start Scan. Either wait a few minutes for the recon queue to decrease or release some of the reserved scan files.
- Bands and lines may be seen in scout images if there is a tube spit or the reference channel is blocked during the acquisition of a scout scan.
- A zero interval axial series, if paused, will add another group equal to the remaining number of images. Make sure to return to the View/Edit screen to delete this group to avoid acquiring additional images.
- Biopsy is not valid for a Thin Twin helical acquisition mode. If Biopsy is selected, the Biopsy window does not open; however a group is added to the series at an inadvertent location. Do not select Biopsy mode with Thin Twin helical scan mode.
- If the message “Can not read Cal Database”, select [End Exam] and re-enter patient information.

SMARTPREP (PURCHASED OPTION)

- Monitor phase scanning will abort if the Monitor ISD is set to less than 3 seconds when the rotation time used is 0.8 sec. Set the Monitor ISD to 3 seconds or greater to avoid this problem.
If this has occurred, do one of the following:
 - Press Monitor Phase to continue monitor scans.
 - Use the graphed data and the clock countdown on the SmartPrep display to determine when to select Scan Phase.
- The film composers will flash and cannot be hidden on the SmartPrep display. This will occur if the Auto Film or Manual Film composer buttons are selected after confirm has been selected for SmartPrep. Don't select the film composer buttons while the SmartPrep display is first coming up or is visible to avoid this.
- If reconstruction fails during SmartPrep, scanning will pause and post a message. Press [Monitor Phase] to resume scanning for the monitor phase.
- SmartPrep may fail to confirm for the Baseline scan. Re –confirm to continue.

PATIENT SCHEDULE

- Close the preferences screen before switching between Patient Schedule and New Patient.
- The Patient Name, Patient ID, Accession Number and Requested Procedure ID can only be edited if “Allow to Edit MWL” is set to yes in the Preference window.

CONNECTPRO (PURCHASED OPTION)

- Bar Codes will fail to be read by the Bar Code reader if the HIS/RIS system the bar code was created on has a different language keyboard than the CT system. For example, if your CT system has a French language keyboard then your HIS/RIS must have a French language keyboard. If it is not possible to have the same language keyboard on each system, then manually enter the Accession or Patient ID number or select the desired patient from the Patient Schedule list to display the patient information on the New Patient screen.
- Patients selected from Modality Work List will be displayed with ^ (carets) to define different DICOM fields. The carets are not displayed on the images.

PERFORMED PROCEDURE STEP (PART OF CONNECT PRO OPTION)

- PPS will post the message “Failed to Start” if the remote server is down.
- Screen Save images created in ExamRx or in the Viewer on the Image Works desktop are not PPS aware in this release of software. Also, images created in Reformat, 3D, and Navigator are not PPS aware. For these image types, “INPR” will be posted in the PPS column in the Browser even though PPS is not enabled.

RECONSTRUCTION

- If images fail to recon the following should be used in order:
 - Select [Recon Management], [Unsuspend Queue]
 - If the image still fails to recon
 - Select [Recon Management], [Restart Queue]
 - If this still fails, Shutdown and reboot the system.
- Lung Algorithm
 - Provides edge enhancement between structures with large density differences, such as calcium and
 - Enhances the contrast of small objects. For best viewing and film quality, select a window width of 1000 to 1500 and a window level of -500 to -600.
 - Increases CT number values at the edge of high contrast objects. When planning to take CT number measurements of vessels or nodules in the lung, please check and compare your results with Standard algorithm images. (ROI and Histogram functions use CT numbers.)
 - The edge enhancement provided by Lung Algorithm may not be appropriate in some clinical cases. Please take individual viewing preferences into account when you choose Lung algorithm.

RETRO RECON

- Two decimal points cannot be entered, i.e. 3.75 for image interval if column edit is used. Edit each scan group (row) if an interval with 2 decimals is required.
- If you delete queued Retro images the image space reserved for those retro recons is not given back. You will not get the space back until you reboot the system.
- Only start Retro Recon when scanning is complete. Do not delete queued retros while scanning is active, scanning could stop.
- If maximum A-P or R-L offset is selected in Retro Recon, one image may fail to reconstruct. Select an offset that is 0.5mm less than the maximum value allowed.
- The start/end locations of a retro reconstruction may be different from Recon 1 and 2 and 3 when a PMR with a thicker slice than Recon 1 is prescribed for the series.
- Retro Recon may hang and core. If the Retro recon screen fails to display or queue retrospective images, retro recon has cored. Reboot the system to recover.
- If an R, L or A, P value is more than half of the Display Field of View in mm. The image annotation will not annotate R and L or A and P on the images. For example if the DFOV is 10cm and the R value is 56mm then the image annotation will show R R. Use a R-L or A-P value which is less than half of the DFOV in mm to avoid this.
- It may be possible to confirm a retro recon prescription with an invalid recon start and end location that results in a recon suspend. To avoid this change image thickness before changing image locations.

RECON MANAGEMENT

- Scout queue entries can be selected in Delete Retro queue entries. Only select Scout queue entries to delete if you are sure you do not need the scout image reconstructed.
- Before saving or restoring scan data make sure the system is idle and no Archive, Network or Filming is active. No other features should be accessed until the save or restore is complete.
- Recon Management may hang while trying to display the menu. Finish the current exam if scanning and shutdown and restart the system to correct the problem.
- If you wish to cancel a Save or Restore scan data, the cancel button will cancel only after the current scan file is saved or restored. A helical scan file contains a large amount of data even though it is one scan file, that file can take up to 30 minutes or more to save. You are unable to scan patients while saving or restoring Scan Data.
Make sure that you have ample time to complete the save or restore before beginning.
- Scan data save or restore is active when the dialog message indicating saving or restoring scan file is displayed.
- If saving scan data, after the save is complete select Restore Scan Files to verify the data is stored.

CT PERFUSION2 ON OC – PURCHASED OPTION

- Saving TIFF files of the graph and images is not supported because there is no floppy drive hardware available on the OC.
- The CT Perfusion Control Panel may be displayed on the current desktop if CT Perfusion is launched from the Browser and the desktop is changed before CT Perfusion screens are completely displayed.

AUTO TRANSFER

- A series may fail to Auto Transfer if [Next Series] then [Create New Series] is selected after a scout scan is prescribed, but not scanned, i.e.; pause and return to view edit is selected.

DISPLAY

- The Exam Rx List Select and Image Works browsers may in some cases not list an exam, series or image. This may be seen after an interruption of power to the system or the database may be updating slowly. The patient does not need to be rescanned.

Verify the following:

- Type the Exam, Series, Image on the accelerator line to display the image,
- Select another exam in the browser this will cause the browser to update,
- The browser may update after a slight delay,
- Reboot the system,
- A second reboot may be needed if a power fail has occurred.

Display cont.

- The browser may update slowly to reflect Screen Save images created in an exam. If the Screen Save images is not listed in the browser wait for them to appear. The time to appear may be up to 15 minutes if the Automatic Database check and recovery is in progress when the Screen Save images are created.
- In general, wait for a display action to complete before entering in another command.
- With the system configured in French or German language, if a comma is entered in explicit magnify, i.e. 2,3 the comma will not be interpreted as a decimal. The magnification factor is applied without the decimal. In the example shown, 2,3 the image would only be magnified 2X not 2.3X. Always enter a period to designate a decimal and have the correct magnification applied.
- Cross Reference lines will post incorrectly if sort by image location is selected. Make sure that sort by number is selected before posting Cross Reference lines.
- Cross Reference lines may not post with some combinations of DFOV and RAS center. RAS centers at the edge of the DFOV and small DFOV's will exhibit the problem. Retro reconstruct the images to a larger DFOV if posting of the cut lines is desired.
- Image annotation may overlap or not be fully shown if a 32 character patient name is entered for 3D, Navigator, Reformat, MID formats, DentaScan images.
- Text and Series pages only display 24 characters for Patient Name.
- Text pages will not show foreign characters entered using the ALT GR key.
- ROI and STD deviation numbers may be reported as zero after a Zoom is applied to an image with a ROI posted. This is due to partial pixel contained in the ROI. Adjust the size of the ROI to include full pixels and recalculate the ROI.
- No text will be entered in the accelerator line if the mouse cursor is over the film composer.
- Some combinations of E/S/I will fail to display an image, if cap locks are on. To select an exam series and image from the accelerator line use the following format when cap locks is on: For example to select Exam 576 Series 2 image 3 type E576 2 3 on the accelerator line.

Selection of an exam series or image may fail a second time after an exam series or image is displayed. If the desired exam series or image does not display then select the desired images from the browser.
- Image location seen on the scanner may differ from that seen on an AW 3.1 system for tilted images. LightSpeed reports the image location as iso-center of the image. AW 3.1 reports image location as the image center for the image. With tilted images there will be a difference in the location numbers proportional to the off center distance times the sin of the tilt angle. The image location can be reported as image center on LightSpeed by using the Report Cursor function.

EXAM RX DISPLAY

- Display may reset itself after a software problem. You may see display stop or not accept your mouse input. Wait for a few seconds and display will automatically reset. You will have to redisplay the images you wish to work on.
- The MIROI pop up may display but will not function. This will occur if an Auto View or Cross reference viewport is in primary focus (blue border) MIROI is not valid in these viewport types. Place a Free or Auto Link viewport in primary focus and reselect the MIROI button to continue.
- Any Report Pixel or MIROI chart that is not screen saved, but is filmed will show the Exam Series and Image as 1000/1/1 on the film composer icon if E/S/I is selected. This is the exam/series/image number the system assigns to this chart display The patients exam number is listed correctly on the filmed image.

Exam RX Display cont.

- Image selection from the accelerator line not function if the primary viewport contains a MIROI plot or Report Pixel chart. Use List Select to display a new image.
- Any accelerator commands entered for Series Binding, Annotation levels for display or filming will be saved as defaults
- ExamRx Display may hang. In some cases switching to the Image Works desktop and back to Exam Rx display may clear the problem. If this doesn't work then perform a system shutdown and reboot the system.
- With paging active (P in the lower left corner of the image) the trackball may adjust window width and window level in other viewports, verify the image has the correct window level before filming the image.
- The right mouse may fail to roam even when the state is set to roam. A system shutdown may be required to clear this problem.
- Display may crash when selecting a 3D object from the list select browser. When a 3D object is selected, 3D should be selected from the browser.
- Images may not be displayed after a switch between AutoView or Image Review layouts. This is due to a restart in the display process. Reselect the desired images for display for the viewport to continue.

IMAGE WORKS DISPLAY

- WW and WL adjustments will only be maintained on individually selected images if the middle mouse button is used. If the accelerator line or Presets are used all images will update regardless of primary or secondary focus of the images.
- Any accelerator commands entered for Series Binding, Annotation levels for display or filming will be applied and not held if the viewer is closed. Use User Preferences in the viewer to save the settings as a default if desired.
- Reference Image may select the wrong scout series, i.e. not the scout series the images were prescribed from if one or more scout series exists. If the wrong scout is selected, place the Reference Image viewport in primary focus and enter the scout series number desired on the accelerator line.
- Window/Level does not display stored value set in the Viewer when +/- series is used. The W/L used is the same as the last series displayed.
- Images created in Add/Sub are displayed with value of WW 4098 and WL 1024.
- If paging is selected while the system is in Compare Mode, the upper left viewport is left blank. Return to the Browser and select the desired series for paging and restart the Viewer.
- Accelerator Command Line entries for E/S/I do not act properly. The first image in the series is always displayed instead of the requested image within the series.
- Cross Reference line numbering along the top of a sagittal scout may be listed in offset manor.
- The new features in User Preferences on the Image Works desktop are not translated.

EDIT PATIENT DATA

- The system may hang if Edit Patient is started while the following operations are in progress:
 - Network Receive
 - Prospective or Retrospective reconstruction is active
 - Archive RestoreTo avoid this, confirm Edit Patient during idle times in network receive, active recon or archive restore.
- Patient age may be seen as zero in Edit Patient Data even though a birth date is entered. Enter the birth date again and patient age will be calculated.
- If duplicate accession numbers exist, Edit Patient Data will not allow the exam to be modified.
- If there is a space at beginning or end of the hospital, Edit Patient Data will fail. Contact your field engineer to update the hospital name correctly in the system reconfiguration.

REFORMAT

- When selecting reformat wait until the % loading dialog is displayed before switching to the Exam RX desktop. This will avoid having to select Load slices to start the reformat model build when returning back to the Image Works desktop.
- The reformat location will not be annotated on images filmed when Printer is selected in Batch Mode. If it is desired to have the reformat location annotated, then use Save in Batch Mode and film the images through Print Series.
- In Batch mode, images reconstructed using “Rotation” or “Oblique” modes will start with number 2 while images created using “First View/Last View” start with 1. In Rotation and Oblique modes, a plan of the reconstructions is screen saved and becomes image 1 for the series.
- Occasionally, the color for the viewports and control panel will change because there is not enough colors available in the resource file to support reformat. To recover, quit reformat and restart the application.
- Window/Level for images are not displayed using the value stored in the image. Use the window/level function keys or selects within Reformat to reset the W/L.

FILMING

- Images or Auto Film control buttons may fail to display in the AutoFilm viewport. This may occur after a switch between desktops. To display the images change to a different AutoView layout, then switch back to your desired layout.
- Don't let the auto film viewport back up with images to be filmed. Start Auto Film as soon as possible to keep caught up on filming.
- Any of the Auto Film control buttons may be activated when the cursor is over the button and the space bar is selected. Place the cursor in the Auto Film viewport only when needed to avoid this problem.
- Auto Film may fail to film. This may occur when the Auto Film viewport selection shows active and the Auto Film status shows paused. Toggle the [Pause] button in the Auto Film viewport and select [Start New Sheet] or [Continue Same Sheet] to restart filming.

Filming cont.

- When scanning multiple one rotation axial groups, images are filmed out of order when Auto Film is selected. Turn Auto Film off for these scan groups.
- Intermittently, Auto Film will automatically print the last page when Add Group is selected.
- Occasionally, the last image will not print. Cancel Auto Film and manually film the image to the Auto Film Composer. If scouts and reference lines were prescribed, they will need to be filmed manually
- Auto Film may fail to display images if recon is having trouble reconstructing the image or if images have failed to install in the database. This may also occur if a large exam has not reconstructed all images and scanning is started on a new exam. The following dialog will be posted if Auto Film cannot display images:
 - Exam 100
 - Series 2
 - Cannot find 10 images
 - Skip Missing images Continue Cancel Film Series

You are presented with 3 choices:

1. Skip Missing images will skip the images that Auto Film cannot find. For example if Auto Film had filmed images 1-10 and you then got the message that 10 images were missing and selected Skip Missing Images, Auto Film would then start filming again with image 21.
 2. Continue will look for the images again, if they are not found then the dialog will be posted again. Before selecting continue verify that recon is active, that the images have been reconstructed, the missing images are not suspended or paused in the Recon Queue and the images are able to be displayed. If the images are suspended or paused in the recon queue, then unsuspend or restart recon to reconstruct the images before selecting Continue.
 3. Cancel Film Series will cancel Auto Film for the series currently being Auto Filmed. The images will then need to be manually filmed.
- If you get this dialog because you have begun an exam while the previous exam still has images to reconstruct, then Select Continue and select Pause Auto Film, resume Auto Film when the Exam Series and Images from the previous exam have reconstructed.
 - When prescribing Mag factor for auto film, scout images if the scout is longer than 500 mm then use a magnification factor less than 1 to display the entire scout.
 - The manual film composer may display when confirm is selected for scanning if the composer had been closed by selecting iconify in the upper right corner of the composer.
 - It's best to run Print Series from only one desktop at a time.
 - It's best to resolve any paused queue entry as soon as possible.
 - Print Series may pause automatically under heavy system load.
 - If Anonymous Patient is selected for an exam where Auto Film is in progress, the system may fail to cross-reference the slices on the scout. If Anonymous Patient is required for series that are still active in AutoFilm, wait until AutoFilm is complete for the series before proceeding with Anonymous Patient.

Filming cont.

- To install a SMPTE pattern select the following list of buttons on the service desktop. Once installed the images will be available for viewing from the Image Works browser or List Select in Exam RX. The SMPTE pattern will be listed as Exam 1000. The patient Name will be listed as SMPTE.
Select Service Desktop
Select [Diagnostics]
Select [Display Processing]
Select [Install SMPTE Image]
- Full Annotation instead of partial annotation will be filmed when using F3 (Film MID) in the Image Works desktop when full annotation is selected.
- The format built in a protocol may be changed when the protocol is used with a message that the format was changed due to an invalid format, even though the format is valid. Verify in Auto Film set up is the format you desire before confirming scan.
- For Multi-Image Display (MID) in AutoFilm, filters, gray scale and orientation selections are applied only to the image in the upper left viewport. If there are an odd number of images in the series, the filters, gray scale and image orientation is applied to the first group only.

ARCHIVE

- The feature status area or the browser may report an Exam has saved even though all images have not been saved to MOD. Scroll through the images in the image window. Highlight the images listed as Archive N and resave them.
- Restoring Exam, Series or images that already exist on the system disk will not post a message that the images are restored or that they already exist. If you have restored images, but get no message that they have been restored, verify that they don't already exist on the system disk.
- When the archive media gets close to being full, the system will always look to see if the exam, series or image can fit on the MOD, if it is desired to have exams sequentially on the MOD, place a new MOD in the drive when first notified the disk may be full.
- If Label is selected in Archive, finish the Label process; don't select cancel on the Label pop-up. Archive will be non-functional; a system reboot will be necessary to recover.
- It's best to resolve any paused queue entry as soon as possible.
- To minimize corruption of MODs, it is extremely important the MOD media be detached and removed from the MOD drive before doing a shutdown or recycling power to the system.
- When restoring from a MOD, it is recommended to write protect the media before placing it in the drive. If the system is unable to read the media, remove the write protect and see if the system is able to perform a recovery of the media to access the data.
- If the MOD cannot be detached or will not eject from the drive, do a shutdown and eject the MOD once the message "Hit any key to start the system" is seen or immediately after the system starts the boot up process. Do not use the screw to eject the disk.
- If you get a timeout message when trying to access the MOD, the MOD drive may no longer be recognized by the system. Perform a shutdown to re-establish communication between the system and the MOD drive.

NETWORK

- LightSpeed images will **not** transfer if Advantage Net protocols are used. Always use DICOM protocol to send LightSpeed images.
- LightSpeed images **cannot** be sent to a HiLight Advantage, HiSpeed Advantage or CT Independent console. These systems do not support DICOM receive.
- Some 3rd party workstations may fail to receive LightSpeed scout images. This is due to the matrix size of the scout image. These stations do not support receiving matrix sizes greater than 512. Some LightSpeed scouts images have matrix sizes greater than 512. If it is desired to have the scout image on the workstation, Screen Save the scout and then transfer the screen save image to the workstation.
- If you query a LightSpeed system from a CT/I or Advantage Windows system in Advantage Net, only Advantage format image exams will be displayed, No LightSpeed exams will be shown. Always query a LightSpeed system using DICOM protocol.
- Lateral scouts displayed on an AW 3.1 system will initially display with zero rotation; they should display with a 270-degree rotation. Rotate the image in the left direction using the rl command on the command line to display the scout in the desired format.
- Images networked to a Advantage Windows 3.1 workstation running software versions prior to 3.1_07 will display a DFOV less than what is displayed on the scanner. This is due to the AW not taking in account the pixels under the focus border of the viewport.
- Images networked to a AW 3.1 will have the following annotation missing:
 - Accession number,
 - Date of Birth,
 - 4i, 2i or 1i recon mode on Axial, Retro Axial or Axial Reformatted images,
 - Series type on Retro images,
 - Table Speed on Helical images,
 - Images are annotated +C when Oral Contrast is used. View the Series Text Page to see if IV contrast was used.
 - CardIQ annotation for BPM and percent of R-R interval and scan type will displayed as Axial.
- The number of images indicated in a series may be incorrect on the remote browser when a query is made from an AW or CT/I station.
- Series types may be listed differently when images are networked to an Advantage Windows 1.2, 2.0 or 3.1 system than that shown on the LightSpeed system
- Exams with Swedish, German or French characters in any Patient Info field will not be transferred to an Advantage Windows.
- When using DICOM protocol, the entered host name must match in spelling and case sensitive, otherwise a connection error message will be displayed.
- It's best to resolve any paused queue entry as soon as possible.

IMAGE MANAGEMENT

- Don't let image space fall below 200 images on a single disk system (130,000-image storage). This will ensure there is room on the system disk to confirm scans and to install reconstructed images.
- Remove images when scanning is idle. This will assure that all images are reconstructed, displayed and installed into the database.

This page left intentionally blank

GE Medical Systems

P.O. BOX 414

MILWAUKEE WISCONSIN 53201-0414 U.S.A.