



GE Healthcare

gehealthcare.com

Technical Publication

**Direction GEHW-SN112017
Revision 2**

NGPDU Gantry Switch Upgrade Procedures

Copyright © 2005 by General Electric Company
All Rights Reserved.



LEGAL NOTES

GE EXCLUSIVE USE ONLY

Restricted To General Electric Medical Systems (GEMS)

This Material is the property of GEMS and contains information that is restricted to use by GEMS Manufacturing, Engineering and Service personnel. Any reproduction in whole or in part is strictly prohibited. Neither this material nor any information contained herein shall be disclosed to anyone outside of GEMS Manufacturing, Engineering and Service personnel.

If you are not a GE Medical Systems Employee, you must exit this material now.

TRADEMARKS

Adobe, the Adobe logo, Acrobat, the Acrobat logo, Exchange, and PostScript are trademarks of Adobe Systems Incorporated or its subsidiaries and may be registered in certain jurisdictions.

Microsoft is a registered trademark and Windows is a trademark of Microsoft Corporation.

All other products and their name brands are trademarks of their respective holders.

COPYRIGHTS

All Material Copyright © 2005 by General Electric Company, Inc. All rights reserved.

Revision History

Revision	Date	Reason for change
1	03/04/2005	Initial release document.
2	04/04/2005	Gantry Switch update kit has been created. Added function test procedures.

NGPDU GANTRY SWITCH UPGRADE PROCEDURES

1 Applied To:

NGPDU, NGPDU-2 and NGPDU-3

2 Problem

As old version NGPDUs have no Gantry Switch, FEs have to open the Gantry cover when they want to turn on/off the Gantry power. To resolve this problem, Gantry Switch is added. This Service Note describes how to install Gantry Switch on old version NGPDU, NGPDU-2, NGPDU-3 which are without this function.

3 Solution


3.1 Overview

Order the Gantry Switch Update Kit Assy (SN: 5136702, refer to [Table 1](#)) and follow the procedures in Section 3.2. [Figure 1](#) gives a pictorial overview of the installation.

Table 1 Gantry Switch Update Kit Assy (SN: 5136702)

ITEM	PART NAME	QTY	PART# (NEW)
1	Front Cover	1	2354578-2
2	Panel Plate	1	2354559-2
3	Real Switch (Gantry Switch)	1	2396924
4	Wire 18# Assy	1	2356917-2
5	Signal Cable	1	2356907-2 (for NGPDU and NGPDU-3)
6	Signal Cable	1	2356915-2 (for NGPDU-2)
7	Clamp	100	U0201AC

Figure 1 NGPDU Gantry Switch Upgrade Overview

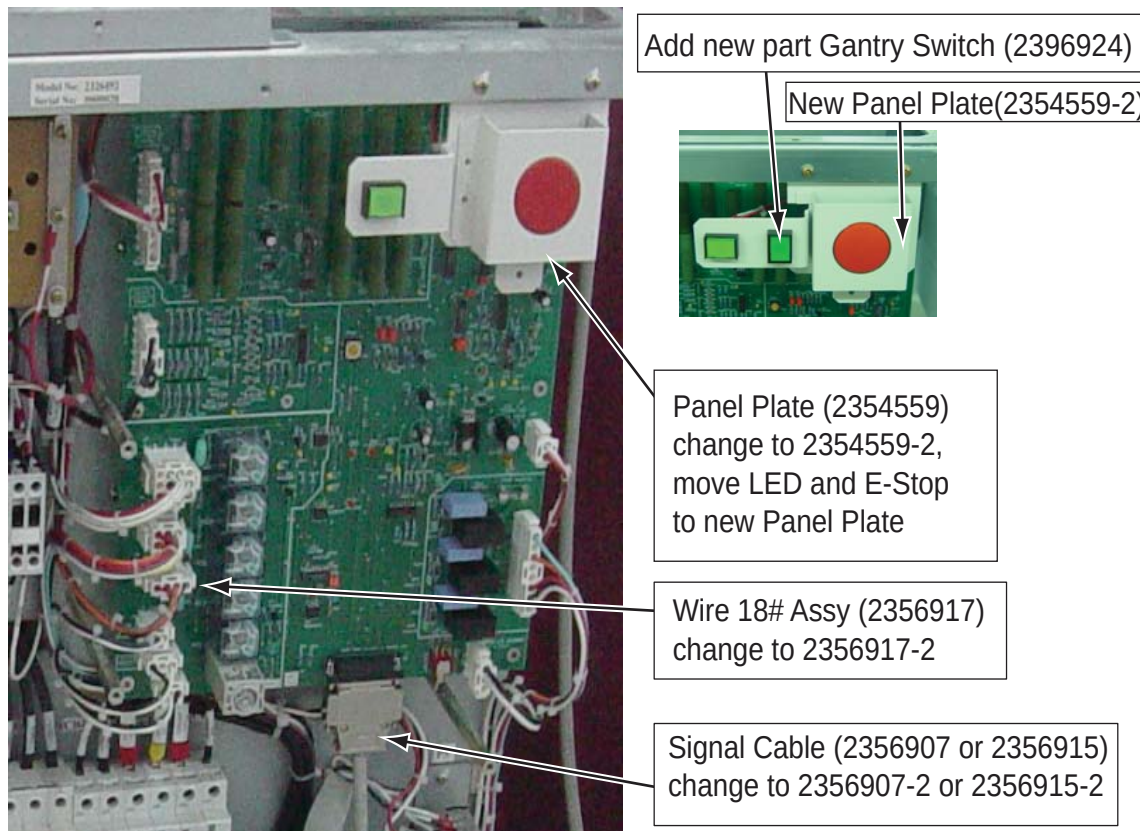

3.2 Procedures



DANGER

POTENTIAL FOR SHOCK.

VOLTAGE MAY BE PRESENT.

TURN OFF AND LOCKOUT THE SYSTEM POWER AT THE MAIN INPUT (POWER DISTRIBUTION BOX).

1.) Preparation

- Turn OFF and lockout the system power at the main input (Power Distribution Box).
- Remove the PDU top cover, front cover and the plastic protection covers.
- Turn off CB1.
- Test the voltage of TS2 (pin1, pin3) with a multimeter and make sure it is less than DC20V (<20VDC).

2.) Gantry Switch Installation

- Remove Panel Plate by unscrewing its two screws. Disconnect Wire 18# and Signal Cable.
- Remove E-stop switch and LED Assy from the original Panel Plate.



Note: No screws need to be removed while removing E-stop switch, but it is very tightly fastened.

- Install original E-stop switch, LED Assy and new part Gantry Switch (2396924) to the new Panel Plate (2354559-2).
- Link Signal Cable (2356907-2 or 2356915-2) to Gantry Switch (2396924), PDU Control

Board and the cable connecting Gantry.

- e.) Connect new Wire 18# Assy (2356917-2) to Emergency Stop Switch, Lamp, Gantry Switch and PDU Control Board.

Figure 2 Connect Cables to Gantry Switch

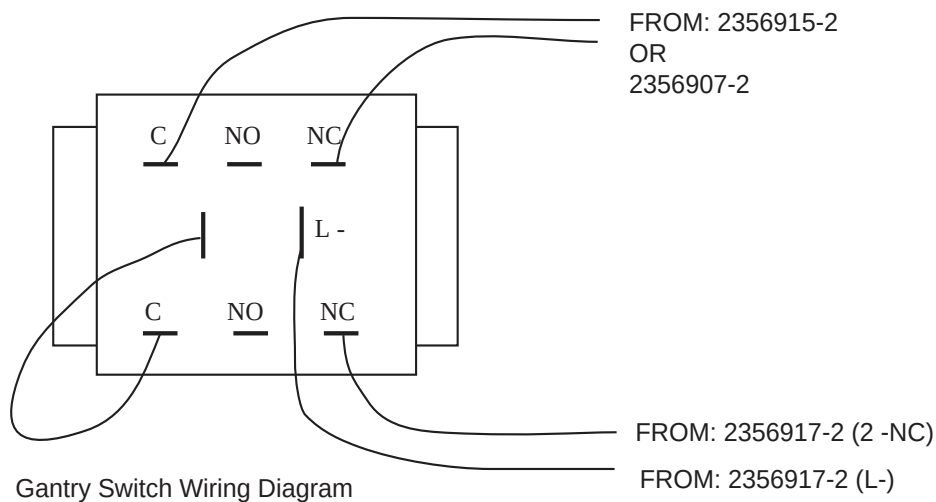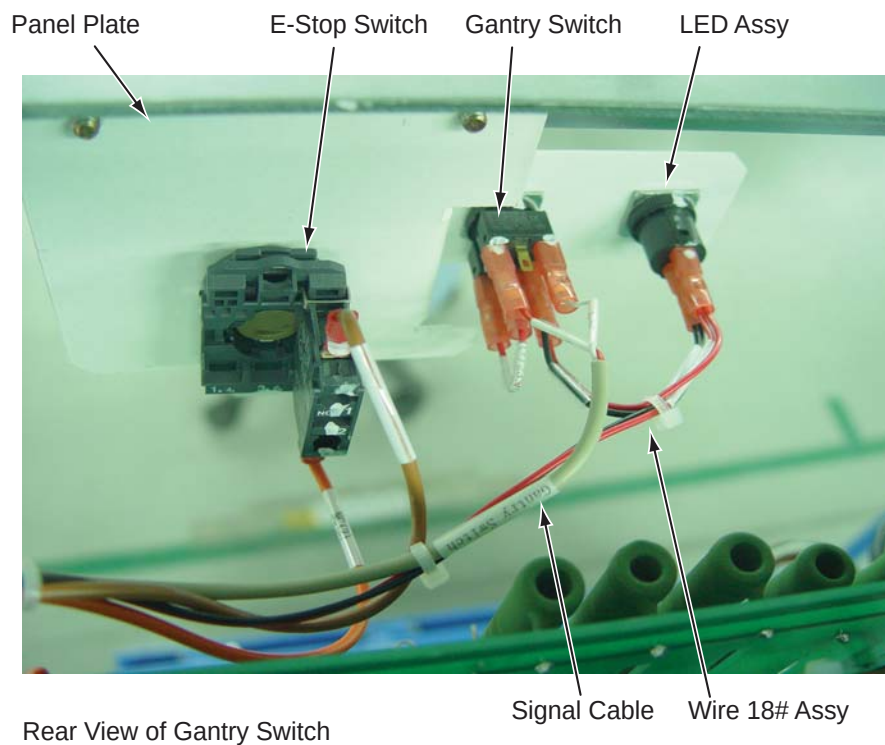

- f.) Install the new Panel Plate (2354559-2) on the PDU.
- g.) Route the cables as shown in [Figure 3](#) and [Figure 4](#). Bundle and fasten the cables with tie-wraps.

Figure 3 Wiring

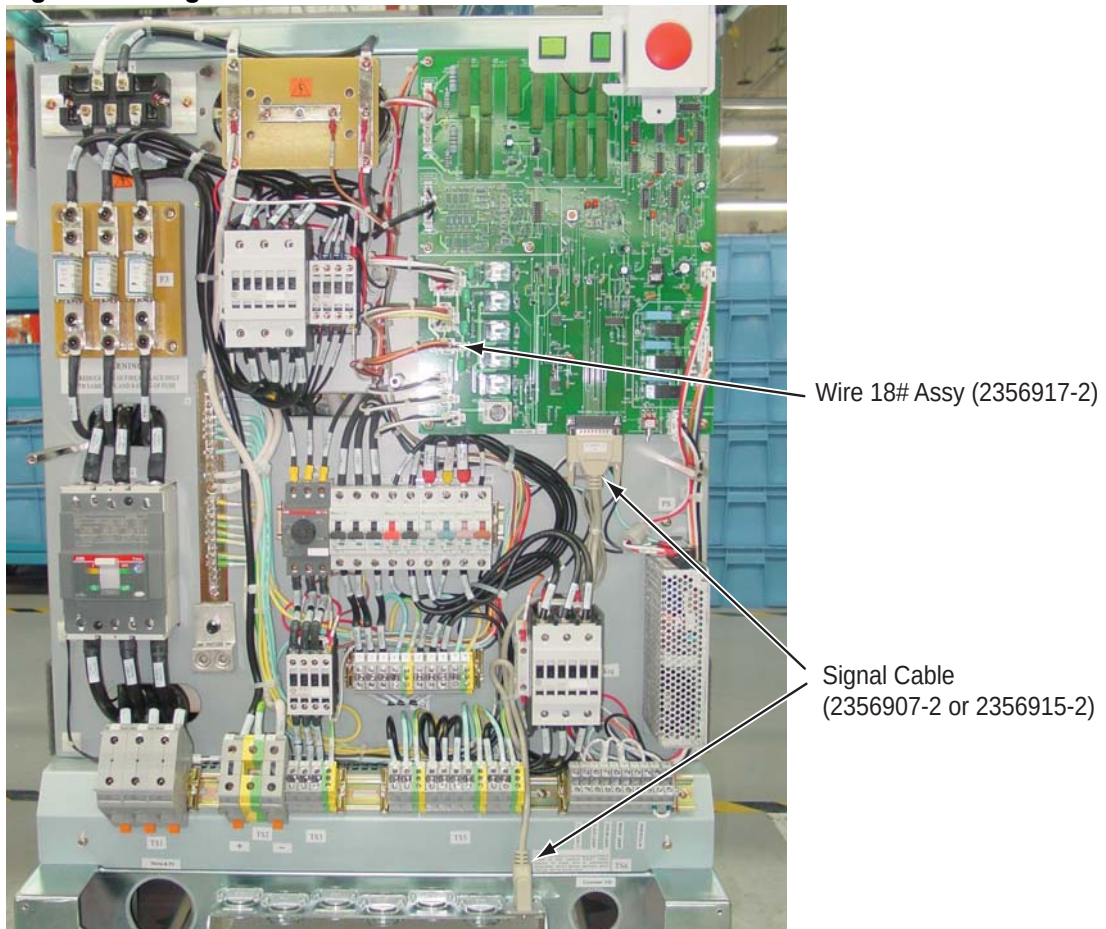


Figure 4 Routing Cables

Top View



-
- h.) Install the protection covers.
 - i.) Turn on CB1.
 - 3.) Function Test
 - a.) Turn on the main power at Power Distribution Box on the wall.
 - b.) Press the Reset button on the Gantry. Verify that OC start up normally and Gantry pass selftest.
 - c.) Press the E-stop Switch. Verify that Table power shut off and Gantry be at emergency stop state.
 - d.) Reset the Gantry again.
 - e.) Press the Gantry Switch. Verify that Gantry Switch LED go out and Gantry power off.
 - f.) Test the PDU output voltage (TS5).
The values of TS5-4 ~TS5-7, TS5-6~TS5-7, TS5-9~TS5-7 should all be about 0V.
 - 4.) Completing Installation
 - a.) Install the new Front Cover (2354578-2) and the top cover.

www.gehealthcare.com



Imagination at work